SERVICE MANUAL 8-CHANNEL SHOULDER COIL M-CAT # M3340CE

NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 1 of 19

WARNING

(EN)

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR, OR PATIENT FROM ELECTRIC SHOCK, OR FROM MECHANICAL OR OTHER HAZARDS.

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NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 2 of 19

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NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 3 of 19

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NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 4 of 19

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NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 5 of 19

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NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 6 of 19

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NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 7 of 19

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NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 8 of 19

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NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 9 of 19

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NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 10 of 19

Revision History

Rev	Date	Author	Primary Reason for Change
A	03/20/2008	G. Wells	Initial revision
B	03/27/2009	G. Wells	Update for MR450/750
C	04/29/2009	G. Wells	Change torque spec requirement note
D	05/06/2009	G. Wells	Add English only warning in 29 languages

Table of Contents

TROUBLE SHOOTING GUIDE.....	13
DIAGNOSTICS	13
SPARE PARTS.....	15
FIELD REPLACEABLE UNITS (FRU).....	16
INSTALLATION PROCEDURES.....	16
COIL MAINTENANCE	19

NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 12 of 19

TROUBLE SHOOTING GUIDE

Follow the table below for possible troubleshooting tips. Otherwise, refer to Instructions for use to contact a service advisor.

Note* There are no service parts in the coil.

Issue	Probable Cause	Resolution
Coil becomes unrecognizable to system	Coil ID issues, bent or broken ID pins	Check Error log and confirm coil ID is read. If not, replace cable and recheck. If still not resolved, replace coil.
Multicoil Bias Failure	DC path broken/open to PIN diodes	Run Diagnostics. *See procedure below.
Coil connector does not fully engage into system interface	Bent pins on connector	Visually check for bent pins. If bent pins are found, replace cable
System housing Red light turns on when coil is plugged into system	Bent or broken contact pins	Visually check for bent pins. If bent pins are found, replace cable
Artifact in image	Magnetic material in scan volume	Remove foreign material from scan room

Table 1. Troubleshooting Guide

DIAGNOSTICS

Gather Load Data – For Excite (14.0 and 15.0) based systems

1. From the Tools menu: Utilities → Service Browser → Diagnostics → Manual Diags → Driver Module
2. Click on Gather Load Data, 5 volts, Enter “5” for file name, Click on “Run”
3. Open new window (Right click) → Service Tools → Command Window
4. Enter into command window: gedit
 - a. Open → Filesystem → usr → g → service → log → 5_mc.txt
5. View open and short circuit results
 - a. Ensure no open circuit values are $\leq -20\text{mA}$
 - b. Ensure no short circuit values are $\geq 200\text{mA}$
6. Run TPS reset after test complete

NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 13 of 19

Multicoil Bias Diagnostics – For Discovery (20.0) based systems

1. From “Common Service” desktop, select Diagnostics → Hardware → Location → Magnet Room → Multicoil Bias Diagnostic.
2. Click “Run”. Diagnostics will indicate pass/fail.

SNR Test [Multi-Coil Quality Assurance (MCQA) Tool]

Description	1.5T HD 8CH Shoulder Coil	3.0T HD 8CH Shoulder Coil
Coil - GE/NeoCoil Part #	5308343-2/NC008901	5271265-2/NC009901
Phantom Holder - GE/NeoCoil Part #	5271265-5/NC008026	5271265-5/NC008026
Phantom - GE Part #	46-265826G6	2359877

Table 2. Tools Required

MCQA Setup and Procedure

**Note – All images below show 3.0T phantoms. 1.5T phantoms contain green solution*

1. Place the base plate on the table and place phantom in the provided slot in the base plate as shown in Figure 1.

Figure 1. Phantom on Phantom Holder Figure 2. Landmark Position

2. Note the Coil Serial Number and place the coil over the phantom as shown in Figure 2.

NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 14 of 19

Note: Ensure the cable rests in the cable guide notch (Figure 1) at the bore-end of the phantom holder and the coil fits snug against notch.

3. Connect the coil to the Hypertronics Connector-A on the Low Profile Carriage Assembly (LPCA) as shown in Figure 3 for HDv systems or as shown in in Figure 4 for HDx/HDx2 systems.

Figure 3. LPCA & Connector for HDv

Figure 4. LPCA & Connector for HDx/HDxt

4. Landmark the phantom at the position shown in Figure 2 and advance to scan.
5. From Common Service Desktop (CSD) select [Image Quality], *Multi-Coil QA Tool* and *Click here to start tool*. Follow applicable service document to run MCQA test. For HDx/HDxt systems instructions can be found in the path:
Troubleshooting → Image Quality → Multi-Coil Quality Assurance Tool.
For Discovery systems instructions can be found in the path:
Troubleshooting → System → Multi-Coil Quality Assurance Tool.

SPARE PARTS

There are no spare parts included with this coil. Field Replaceable Units may be purchased.

NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 15 of 19

FIELD REPLACEABLE UNITS (FRU)

Description	GEHC Part #	NeoCoil Part #
1.5T HD 8Ch Shoulder Coil Assembly with Cable	5308343-2	NC008930
3.0T HD 8Ch Shoulder Coil Assembly with Cable	5271265-2	NC009930
HD 8Ch Shoulder Coil Cable Assembly	5271265-3	NC008931
8Ch Shoulder Coil Phantom Holder Assembly	5271265-5	NC008932

Table 2. FRU listings

INSTALLATION PROCEDURES

8-channel Shoulder coil system cable replacement (GEHC Part# 5271265-3)

Tools required:

Screw driver (Supplier recommended Torque screw driver set to 8in-lbs)

Anti-static station

1. Remove screws from preamp housing.
2. Carefully open housing halves and set aside cover.
3. Remove four (4) cable clamp screws from housing (See Figure 1)
4. Remove screws from four corner locations on the large PCB (See Figure 2).
5. Remove all eight cable connectors (See Figure 3)
 - a. Pull cable connector straight from socket to ensure no damage is done to the socket center pin. Check center pin for damage.
6. Remove DC connector (See Figure 4)
 - a. Use a flat screw driver to gently pry each side until loose
7. Discard cable.
8. Following color code in table 3, attach new cable connectors (See Figure 3) and DC connector (See Figure 4).
9. Set PCB assembly into housing and replace screws (See Figure 5 – Supplier recommended 8in/lbs torque – DO NOT OVER TIGHTEN).
 - a. Pay close attention positioning of cables as stated in Figure 5
10. Replace cable clamp screws. (Supplier recommended 8in/lbs torque – DO NOT OVER TIGHTEN).

NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 16 of 19

11. Carefully replace cover ensuring no cables or wires are pinched by the housing.
12. Replace housing screws. (Supplier recommended 8in/lbs torque – DO NOT OVER TIGHTEN).

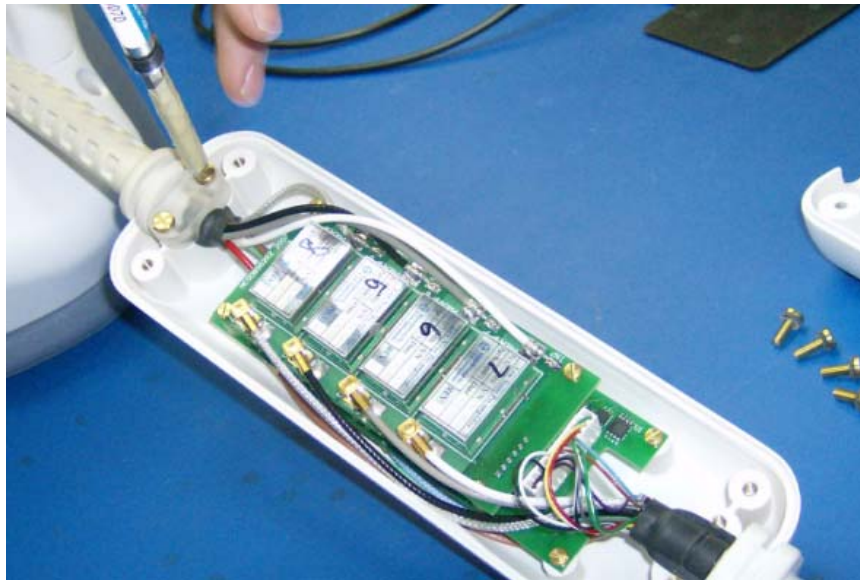
Figure 1. Cable clamp screws

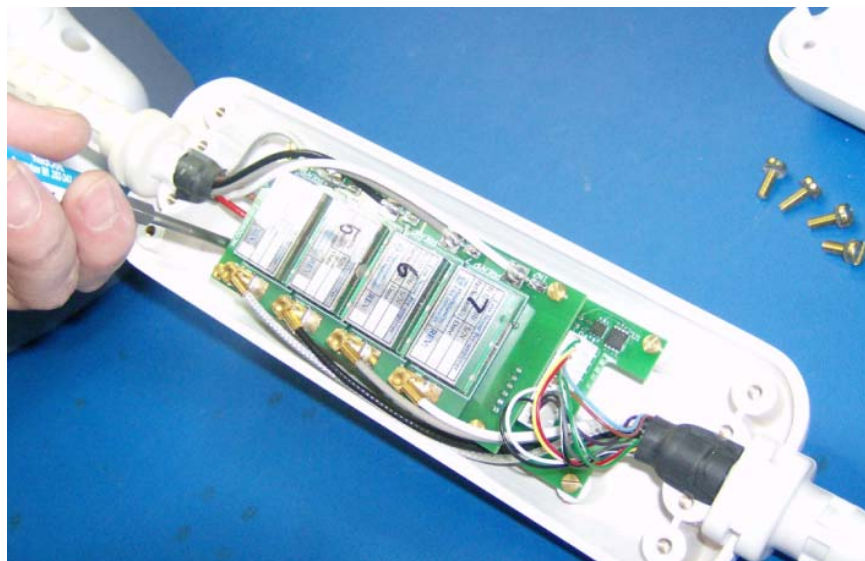
Figure 2. PCB screws

NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 17 of 19

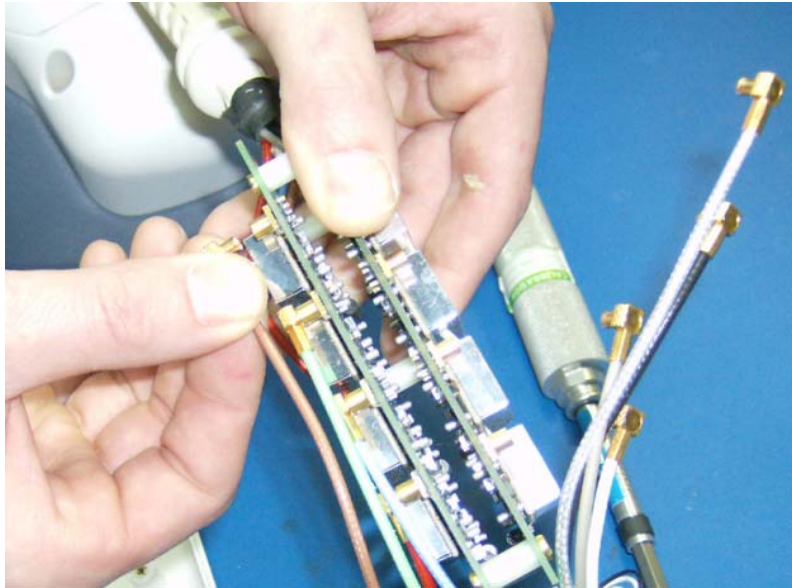
Figure 3. Cable connector removal/replacement

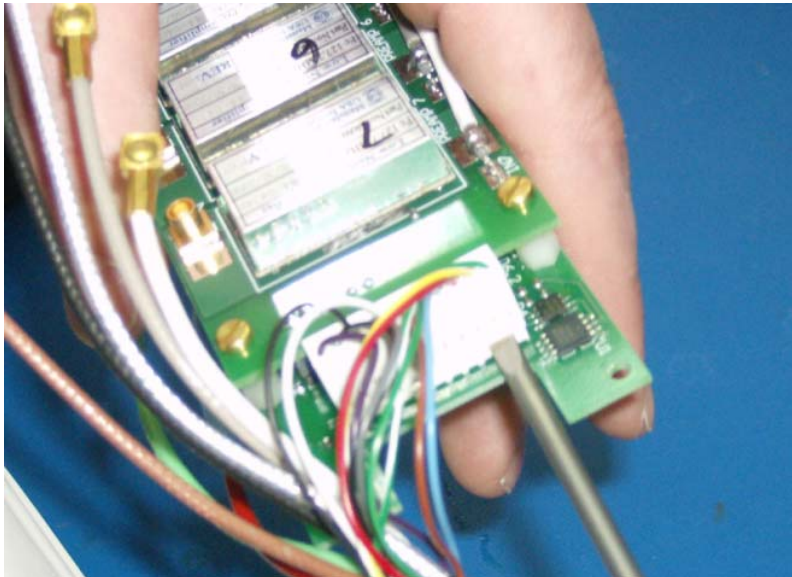
Figure 4. DC connector

NC008SVM Rev D	Service Manual 8-ch Shoulder Coil
© 2008 NeoCoil, LLC NeoCoil Proprietary and Confidential Information	Page 18 of 19

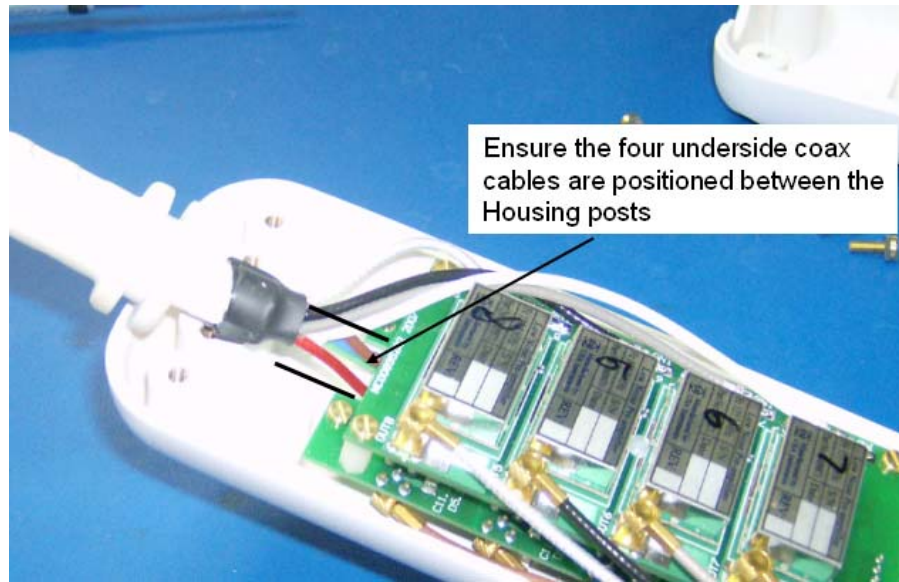
Figure 5. Coil cable placement

Channel	Coax Cable
1	BROWN
2	RED
3	GREEN
4	BLUE
5	BLACK
6	GRAY
7	WHITE
8	CLEAR

Table 3. System Cable Color Code

COIL MAINTENANCE

See Instructions for Use manual