GE Healthcare

Discovery MR450 Mechanical Dock Upgrade

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File a report with

- Call 1–800–548–3366 and use option 8.
- Contact your local service coordinator for more information on this process.

Rev. 08/15/2003

© 2009 General Electric Company
GE Medical Systems, a General Electric Company, going to market as GE Healthcare
3000 N. Grandview Boulevard
Waukesha, Wisconsin 53188
USA

www.gehealthcare.com

TABLE OF CONTENTS

<u>SECTION</u>	<u>PAGE</u>
TABLE OF CONTENTS	1
REVISION HISTORY	3
LANGUAGE POLICY	5
SECTION 1 – GETTING STARTED	
1-1 INTRODUCTION	15
1-2 REQUIRED TOOLS	15
1-3 DAMAGE IN TRANSPORTATION	15
1-4 PRODUCT LOCATOR	16
SECTION 2 – INSTALLATION PROCEDURE	
2-1 SAFETY	17
2-2 DOCK REMOVAL	17
2-3 INSTALLATION OF NEW 32-CHANNEL DOCK	20
2-4 ALTERNATE CABLE INSTALLATION	23
2-5 POWER UP	23
SECTION 3 – COMPLETION CHECKLIST	
3-1 MECHANICAL INSTALLATION COMPLETION	24
3-2 GE FIELD ENGINEER RESPONSIBILITIES	24

This page left blank intentionally

REVISION HISTORY

REV	DATE	PRIMARY REASON FOR CHANGE
-----	------	---------------------------

1	21 Apr 2009	Preliminary release for verification.
---	-------------	---------------------------------------

This page left blank intentionally

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WARNING

(EN)

- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR, OR PATIENT FROM ELECTRIC SHOCK, OR FROM MECHANICAL OR OTHER HAZARDS.

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This page left blank intentionally.

SECTION 1 – GETTING STARTED

1-1 INTRODUCTION

This manual provides instructions to install the Mechanical Dock (or 32-Channel Dock) on the MR450 1.5T System.

This Direction covers the following catalog:

M7000EM – Upgrade for Configuration A to B

1-2 REQUIRED TOOLS

FERROUS MATERIAL HAZARD! IF MAGNET IS ENERGIZED, TOOLS AND PARTS REQUIRED FOR THIS INSTALLATION THAT CONTAIN FERROUS MATERIAL AND WILL BE STRONGLY ATTRACTED TO MAGNET AND MAY BECOME DANGEROUS PROJECTILES. KEEP ALL FERROUS TOOLS AT LEAST 10 FEET AWAY FROM THE MAGNET.

<u>ITEM</u>	<u>QTY</u>
Non-ferrous adjustable wrench	1
Non-ferrous flathead screwdriver	1

1-3 DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation “**damage in shipment**” written on **all** copies of the freight or express bill **before** delivery is accepted or “signed for” by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

File a report with

- Call 1-800-548-3366 and use option 8.
- Contact your local service coordinator for more information on this process.

1-4 PRODUCT LOCATOR

The Global Install Base Database (also known as Product Locator System) tracks shipment, trans-shipment, and field location of the serialized models. There are now two methods for submitting Product Locator information.

1. At this time, for **U.S. ONLY**, the preferred method of submitting information is the *FE Site Verification Web Site* at:
<http://gein2.med.ge.com/gib>

The FE Site Verification consists of three components that are available on the web from the main menu. They are:

- Install/deinstall product locator model and serial numbers
- Add/modify ship to address information
- Update CARES FE data for primary/secondary FEs

To obtain a copy of the FE training tutorial for using this website, a downloadable copy is available at the following:
Product Locator Support Central Page – <http://supportcentral.ge.com/15563>

2. One “Shipping Card” is filled out and submitted when shipped (extra cards are supplied for trans-shipments between storage and distribution points), and the “Installation Card” and extra shipment cards are attached.

Verify that serial and model number on each rating plate matches installation card numbers before removing installation card. Note that there may be one or more shipment cards and bar code labels with the installation card. These shipment cards are used to trace the transfer of serialized units between various inventory storage and distribution points until the product reaches its final installation destination. Process just the installation card and discard any extra shipment cards and labels.

SECTION 2 – INSTALLATION PROCEDURE

2-1 SAFETY

⚠ DANGER

FATAL ELECTRIC SHOCK HAZARD!!

TO PREVENT FATAL ELECTRIC SHOCK, DISCONNECT POWER FROM THE PDU BEFORE YOU PERFORM THE FOLLOWING PROCEDURES.

PERFORM LOCKOUT / TAGOUT PROCEDURE PER GE OSHA LOCKOUT / TAGOUT REQUIREMENTS LISTED IN THE FIELD SERVICE EHS MANUAL & TRAINING 2000 P/N 2135387-200.

⚠ WARNING

POSSIBLE PERSONAL INJURY!

THIS PROCEDURE REQUIRES REMOVAL OF FERROUS MATERIAL IN CLOSE PROXIMITY TO THE MAGNET.



THEREFORE, TWO FIELD ENGINEERS ARE REQUIRED AT ALL TIMES.

2-2 DOCK REMOVAL

1. Perform Lock Out Tag Out.
2. Undock the patient table and remove it from the scan room.


2-2-1 DISCONNECT CABLES

1. Disconnect the three cables routed from the Dock that are connected to the bracket on the Left front leg.
2. If required, unfasten the magnet enclosure cover from the bracket to aid in sliding out the cables.

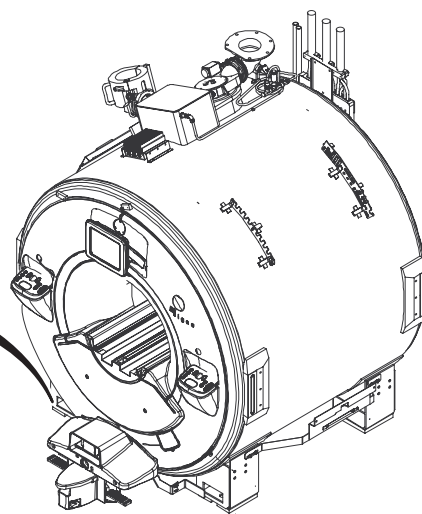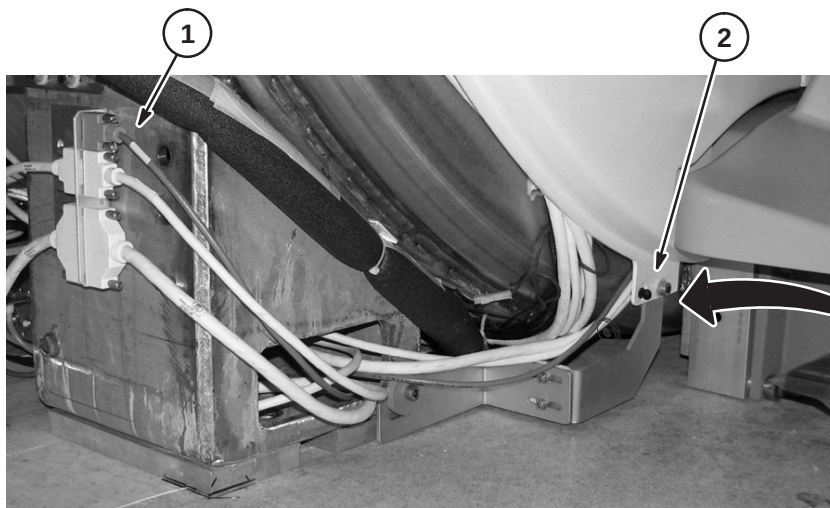
2-2-2 UNFASTEN AND REMOVE DOCK

1. Scribe a line around the base of the dock on the floor to ensure the replacement dock will be located in the same position.

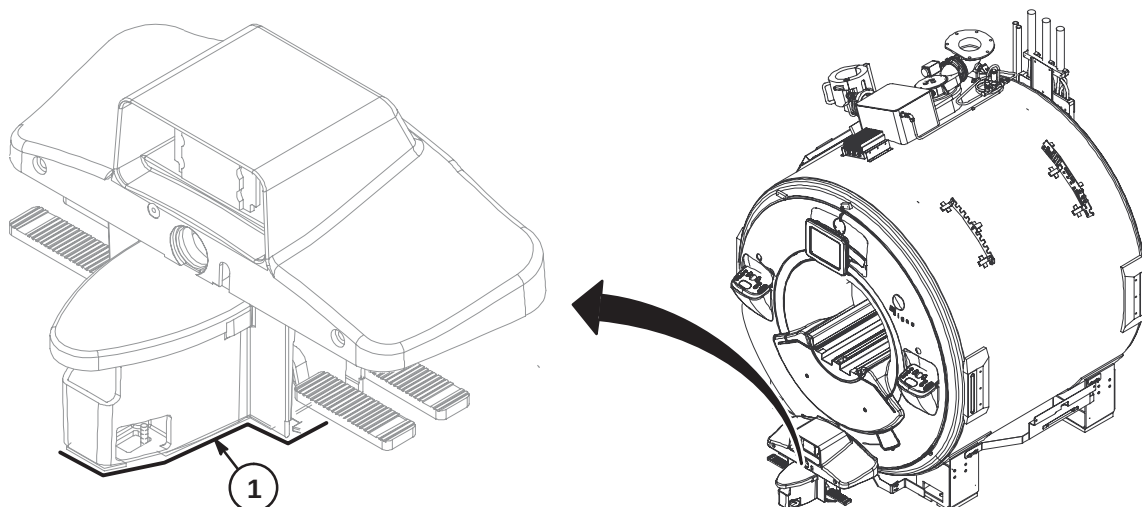2. Loosen the nuts attaching the dock assembly at the front and at the right and left brackets.
3. Make sure the same quantity of washers remain on mounting studs between Dock and mounting brackets when Dock is re-attached.

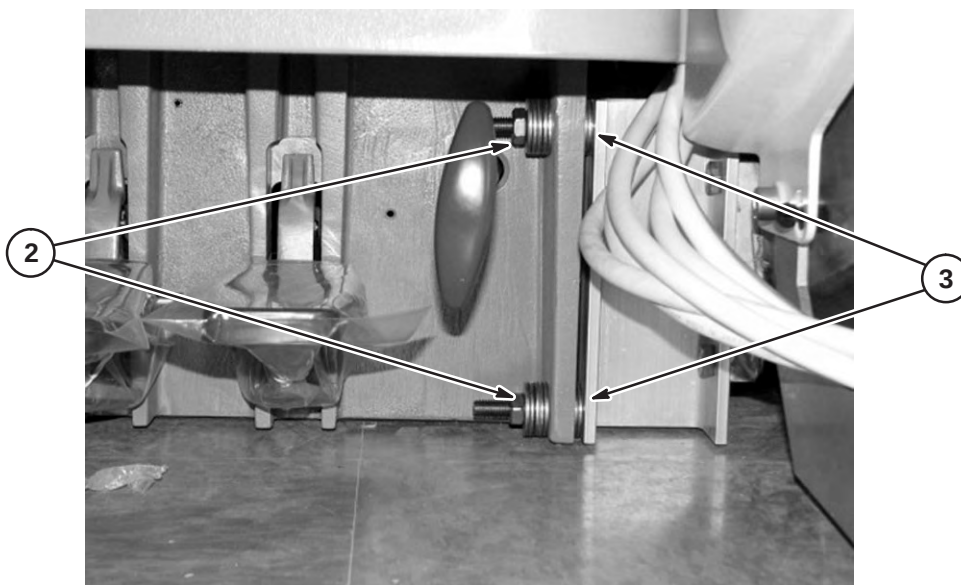
WARNING

POSSIBLE PERSONAL INJURY!
THE DOCK CONTAINS FERROUS MATERIALS.
FROM THIS POINT OF THE PROCEDURE FORWARD, DOWNWARD PRESSURE MUST BE CONTINUOUSLY APPLIED TO THE DOCK SO AS TO PREVENT IT FROM BEING PULLED INTO THE MAGNET BORE.

2-2-2 UNFASTEN AND REMOVE DOCK (Continued)

4. Applying pressure to the top of the dock, completely remove all dock attachment nuts and washers.
5. Use two FEs to pull the dock away from the magnet until the stud in the front is hitting the inner frame.


NOTE: Do not remove washers that are located behind the Dock, as they are needed to proper Dock alignment.

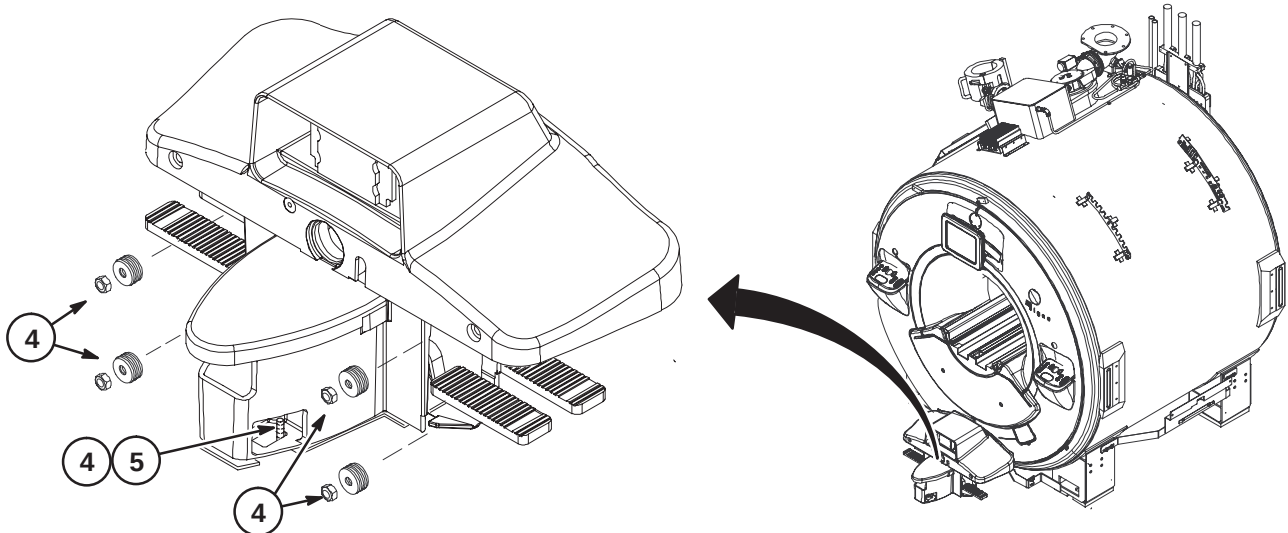**⚠ WARNING****POSSIBLE PERSONAL INJURY!**

TO OPPOSE ANY PULL FROM THE MAGNET, FIELD ENGINEERS MUST APPLY CONTINUOUS DOWNWARD PRESSURE TO THE TOP OF THE DOCK EVEN WHEN LIFTING THE DOCK IN THIS STEP.

DO NOT PLACE ANY PART OF THE HUMAN BODY IN THE PATH BETWEEN THE DOCK AND THE MAGNET.

6. Use two FEs to lift the dock up JUST ENOUGH to clear the top of the stud in the floor.

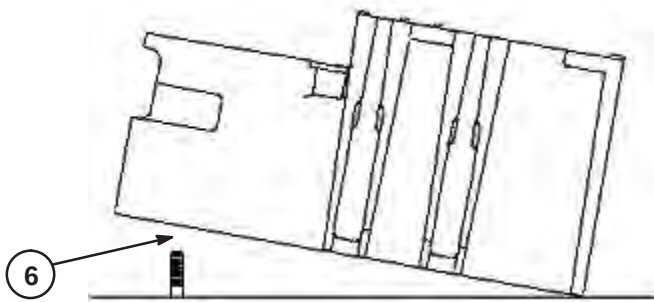
NOTE: It should not be extremely difficult for two FEs to pull the dock away from the magnet once all fastening hardware is removed. If you do encounter difficulty, check that the cables are not caught against the front of the magnet.

7. With two field engineers, one holding each side of the dock, slide the dock away from the magnet front to where the foot of the patient table would normally be located.
8. When a safe distance away from the magnetic field (at least as far as the foot end of the patient table), two field engineers may slowly lift the dock from the floor and exit the room.


2-3 INSTALLATION OF NEW 32-CHANNEL DOCK

2-3-1 INSTALL DOCK


1. With two field engineers, one holding each side of the dock, enter the magnet room and place the dock on the floor where the foot of the patient table would normally be located.
2. Position Rubber Pad under Dock.
3. With two field engineers, one applying vertical downward pressure to each side of the dock, slide the dock toward the magnet front.
4. Make sure cables are routed to each side of mounting brackets while sliding Dock into position.

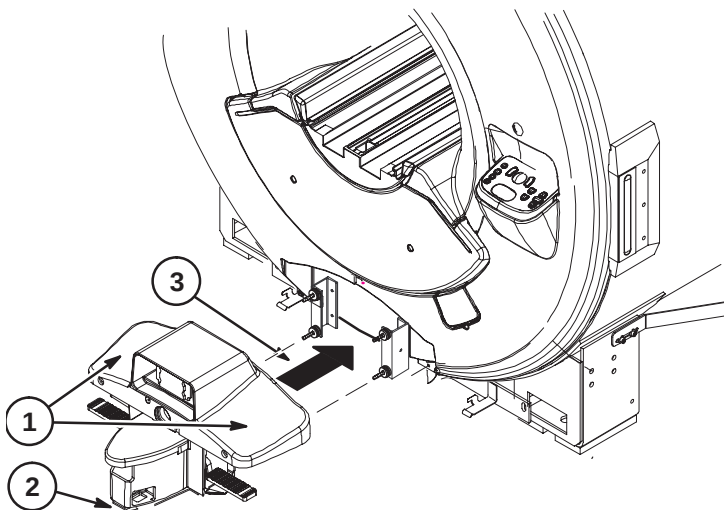5. While maintaining pressure to the top of the dock, position the dock against the base of the magnet to the left of the floor stud at the front of the dock.

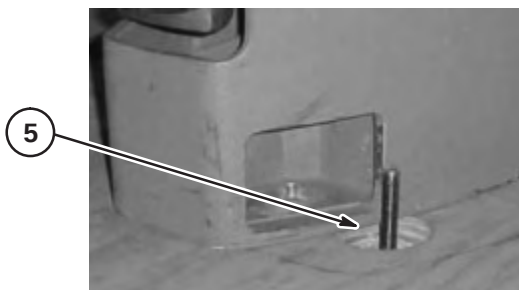
6. Lift the dock just enough to clear the top of the stud in the floor.
7. While rotating the dock to the right, align the hole at the base of the dock with the stud and place the dock down into position.

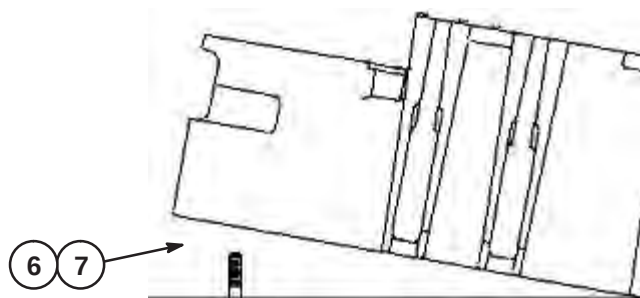
2-3-1 INSTALL DOCK (Continued)

8. Position replacement Dock on scribed line.
9. Fasten Dock Asm to floor mounted stud with previously removed Dock Clamp Plate, 10.5mm ID Washer , and 10mm Hex Nut. Apply Loctite 242 as required.

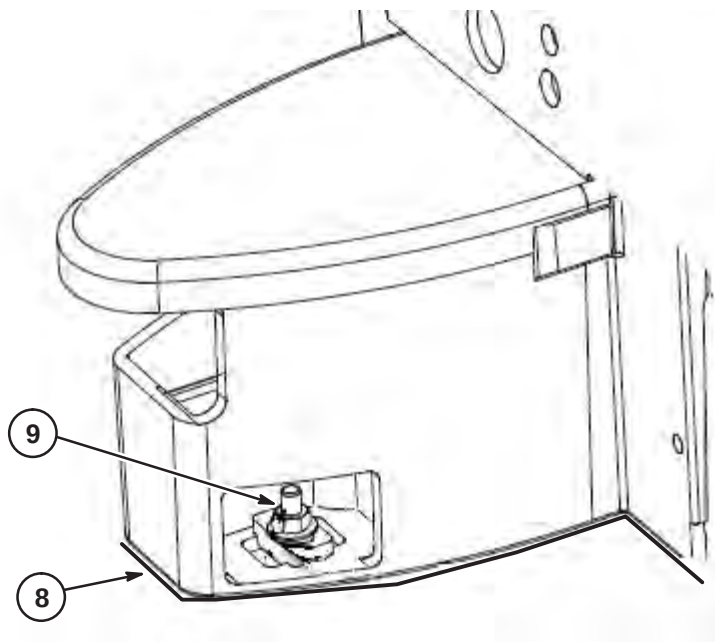10. Secure Dock to mounting studs using previously removed nuts and washers.

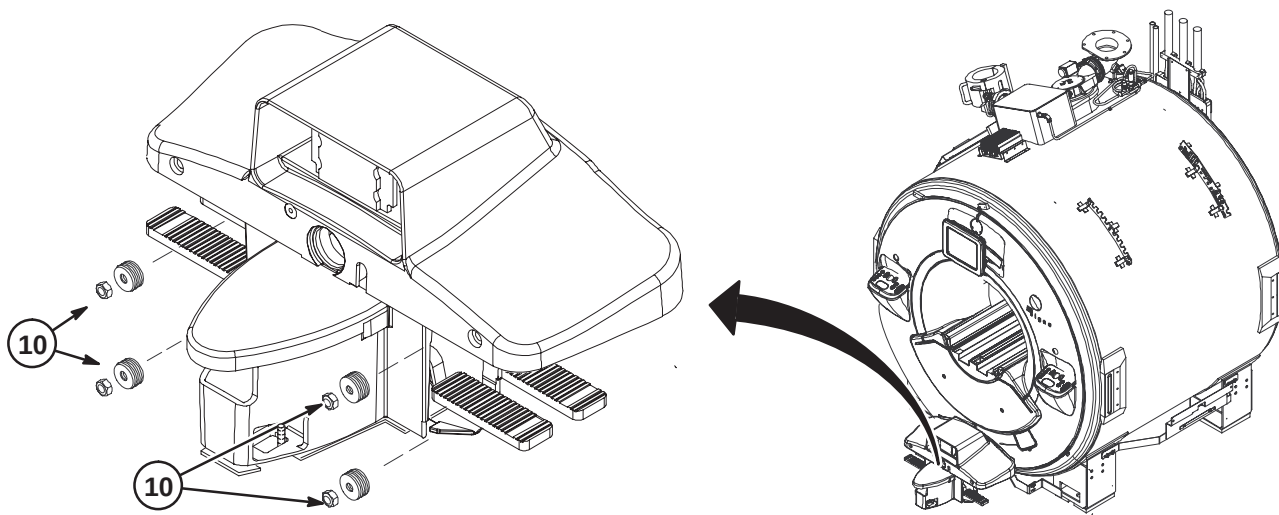
2-3-2 DOCK CABLES – LEFT SIDE

1. Route and connect cables from dock as shown.
2. Insert and loop excess cable thru hole at front of magnet leg.
3. Loosen bottom of cover and slide cables behind cover. This is required to install front Trim Covers.


**2-3-3 DOCK CABLES – RIGHT SIDE**

1. Route cables from dock as shown. Cables can be routed using any of the following three ways:
 - Thru the hole in the magnet (shown below)
 - Between the leg and the magnet
2. Loosen bottom of cover and slide cables behind cover. This is required to install front Trim Covers.


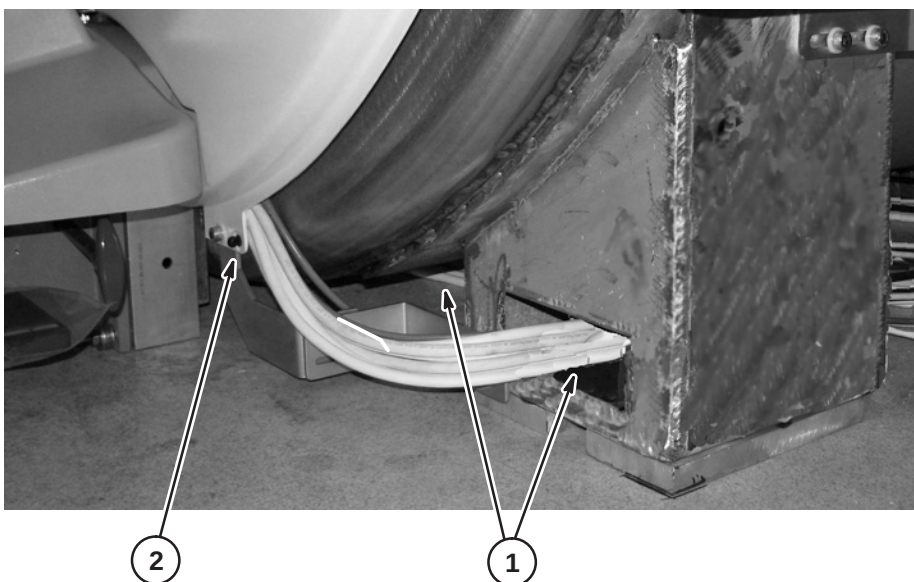
2-3-3 DOCK CABLES – RIGHT SIDE (Continued)

3. Route cables under magnet to Rear Pedestal. Dress and neatly store excess cable under magnet.

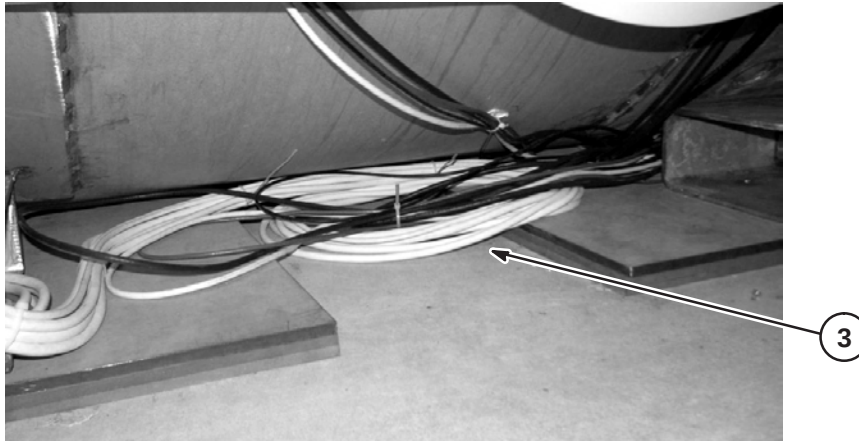4. To right side of Rear Pedestal – route and connect Runs M1330, M1331, M1332, and M1333 to PED-RBOX as marked.
5. To left side of Rear Pedestal – route and connect Runs M1371 and M1372 to SLOT 11, J3 & J4 as marked.
6. To left side of Rear Pedestal – route and connect Run M3373 to J9 on PED-RBOX as marked.

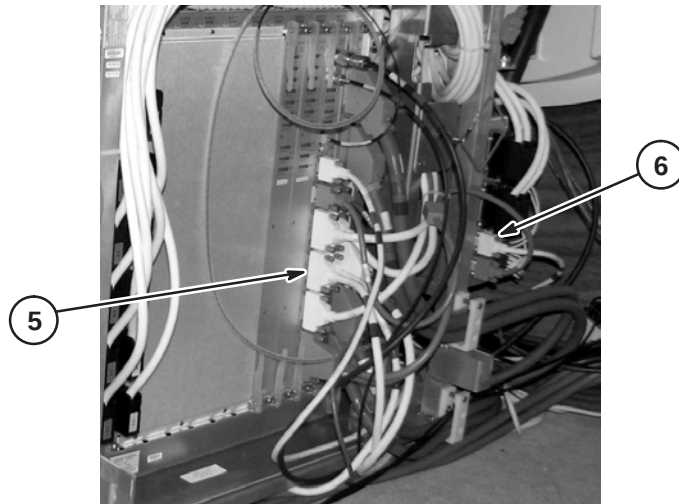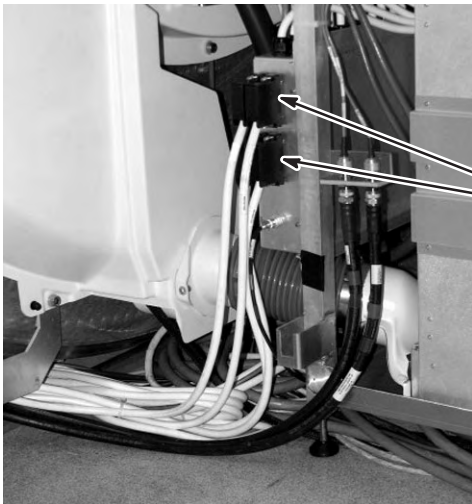2-4 ALTERNATE CABLE INSTALLATION

If required, an alternate set of cables (5192500) is provided for routing and connecting from the cable bracket on the side of the right front magnet leg to the Rear Pedestal.

2-5 POWER UP

1. Remove the power lock-out tag-out and restore power.
2. Align and dock the table.

Note: If adjustment is required for alignment of patient table, refer to calibration and setup procedures on Service Methods DVD until proper gap is achieved.

SECTION 3 – COMPLETION CHECKLIST

3-1 MECHANICAL INSTALLATION COMPLETION

The following steps need to be reviewed and completed by GE Field Engineer, especially if the system was installed by non-GE personel. Check off boxes at the front of each item below when complete.

- ☐ Check procedure to make sure all steps are complete
- ☐ Resolve shipment shortages
- ☐ Make sure all rating plates are installed
- ☐ Recheck all wiring connections using cable map
- ☐ Correct for wiring errors if necessary

3-2 GE FIELD ENGINEER RESPONSIBILITIES

1. Record and enter applicable data into applicable site configuration files and records.
2. Complete Product Locator information for all installed serialized components, new or updated, via either:
 - **(U.S. Only)** *FE Site Verification Web Site*
or
 - Process and return product locator installation cards for all serialized components to:
Product Locator File, P.O. Box 414, W-523, Milwaukee, WI 53201-0414

Refer to Section 1-4 for details on submitting Product Locator information.

Note

Failure to fill out and return Product Locator Cards may result in failure of your site to receive future FMI's.

3. Store the delivered site's set of service tools and spares kit in service cabinet at site.
4. Leave the site's set of Manuals and CD-ROM at site. Install Service and Operator documentation on the host computer. Organize and set up reference cabinet for the Manuals.
5. Locate any Material Safety Data Sheets (MSDS).. They must be retained on site. Customer is to be informed that material with MSDS was brought into site and customer should know/decide where MSDS should be retained at site.