



3 PRE-INSTALLATION REQUIREMENTS

3-1 EQUIPMENT PLACEMENT

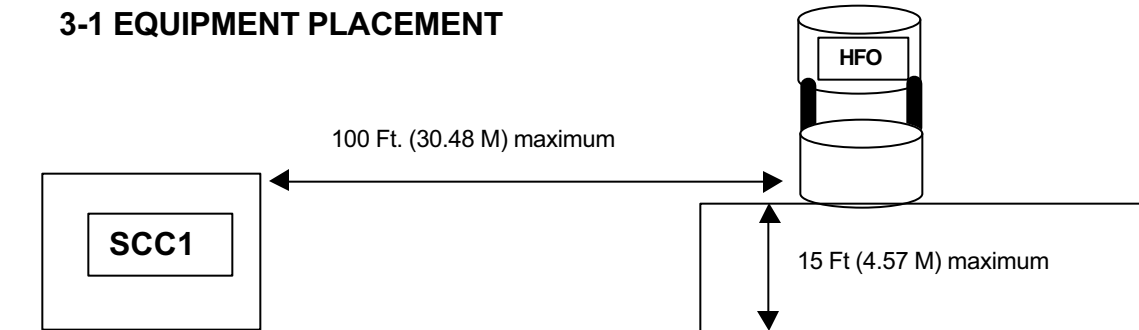

The maximum pipe length from the SCC1 to the magnet is 15 feet (4.57 m) vertical and 100 feet (30.48 m) total.

Illustration 1 **PLACEMENT OF SCC1 TO MAGNET**

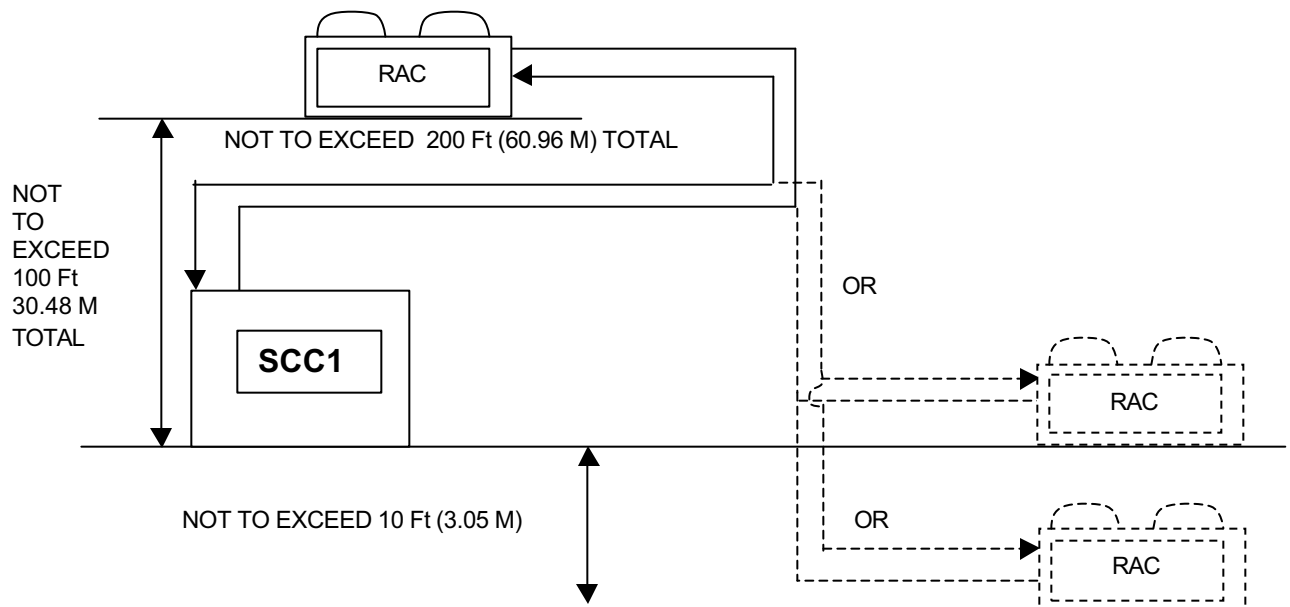

The SCC1 and the RAC pipe length should not be separated by a distance greater than 200 equivalent feet (60.96 M).

Illustration 2 **PLACEMENT OF SCC1 TO RAC**

3-2 SITE ENVIRONMENT**3-2-1 Operating Temperature And Humidity Specifications**

Item	Requirement
Ambient Temperature	41°F (5°C) to 104°F (40°C)
Humidity	0% to 100% RH
Room Cooling Requirements	12,000 Btu/hr
Altitude	984 feet (300 m) below sea level 11,060 feet (3,597 m) above sea level
SCC1 and SCC2 Noise Level Maximum	75 dBA at 1 meter in all directions
Magnetic Field	30 Gauss

Table 5 **Operational Conditions****3-2-2 Power Requirements**

Item	Requirement
Wire - NOTE:	4 AWG (25 mm ²) minimum 3 conductor w/ground/ 2 AWG (35mm ²) maximum
Power Consumption SCC1/SCC2/RAC Note	39 KVA (maximum)
Voltage Requirements	480 VAC ± 10% 60 Hz ± 3 Hz 380 VAC ± 10% 50 Hz ± 3 Hz
Circuit Size	90 AMP (maximum)
NOTE: The SCC1 will be hardwired to the MDP (Main Distribution Panel) with a four conductor (4 AWG / 25 mm ²) power cable supplied by the installer. The Shield Cooler and Cryo Cooler compressors are part of the power consumption.	

Table 6 **Power Requirements**

3-3 EQUIPMENT CLEARANCE REQUIREMENTS

It is necessary to maintain proper distances at the top, front and back of the SCC1 cabinet in order to have a good service area and cooling clearance. Consideration has to be made to keep the area above the refrigeration piping clear if the piping needs to penetrate the ceiling/roof area. There is no requirement for clearance at the side of the cabinet, but any allowances would aid in the service of the SCC1 cabinet.

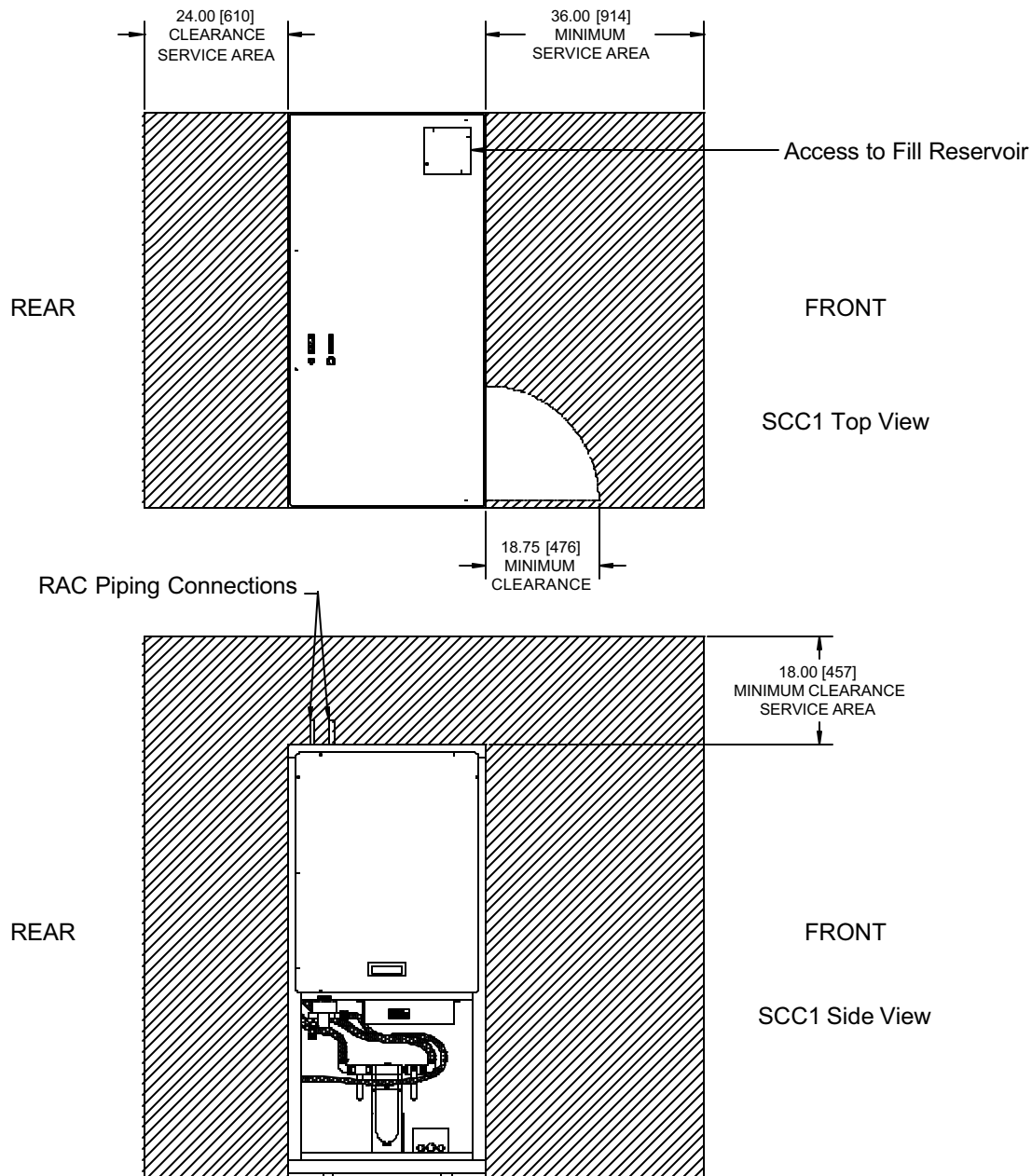


Illustration 3 **NECESSARY CLEARANCES FOR SCC1 ([] = cm)**

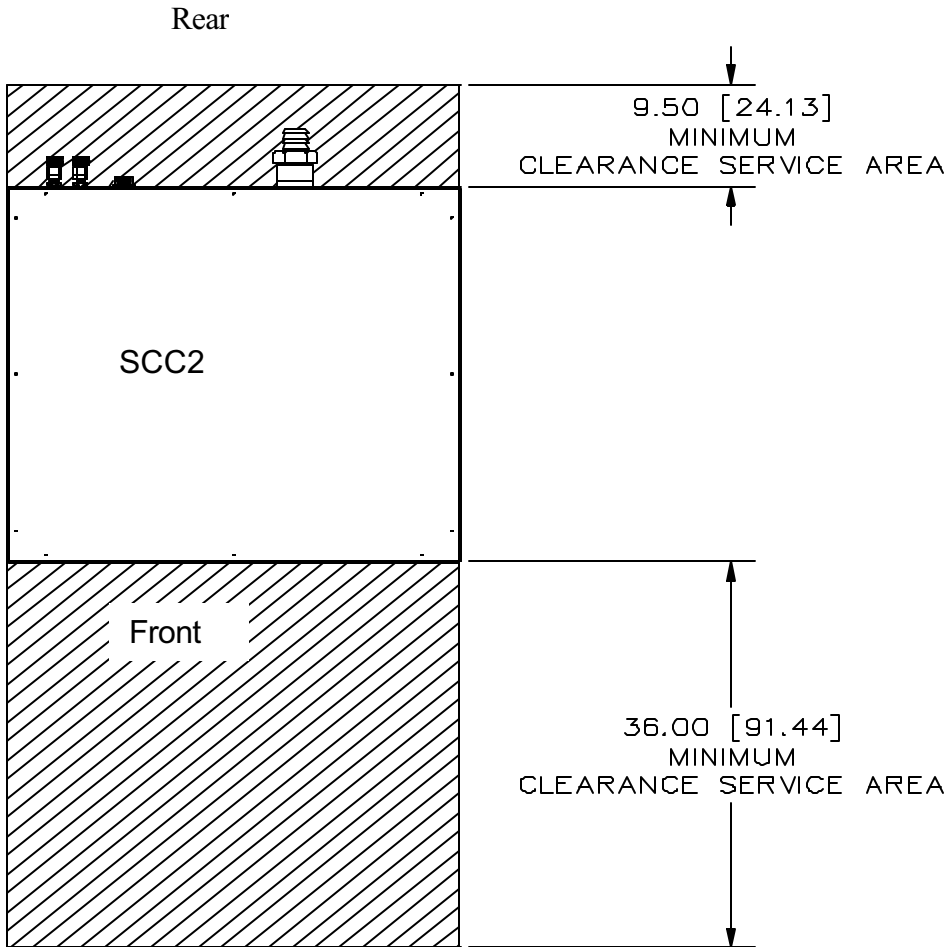


Illustration 4 **NECESSARY CLEARANCES FOR SCC2 ([] = cm)**

Minimum clearance service area is required to allow cabinet to be rolled forward for increased service access at rear.

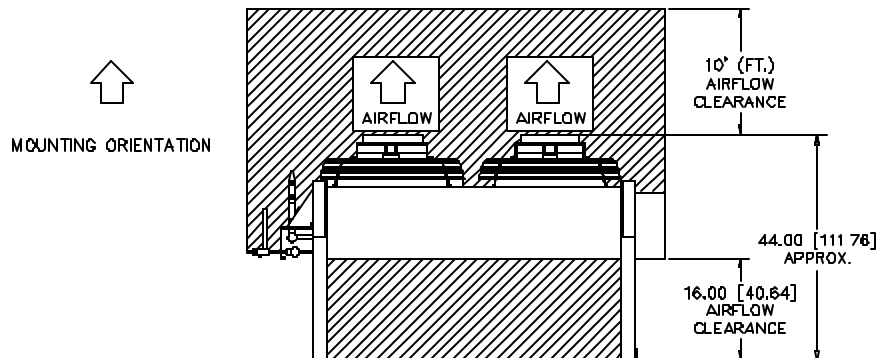
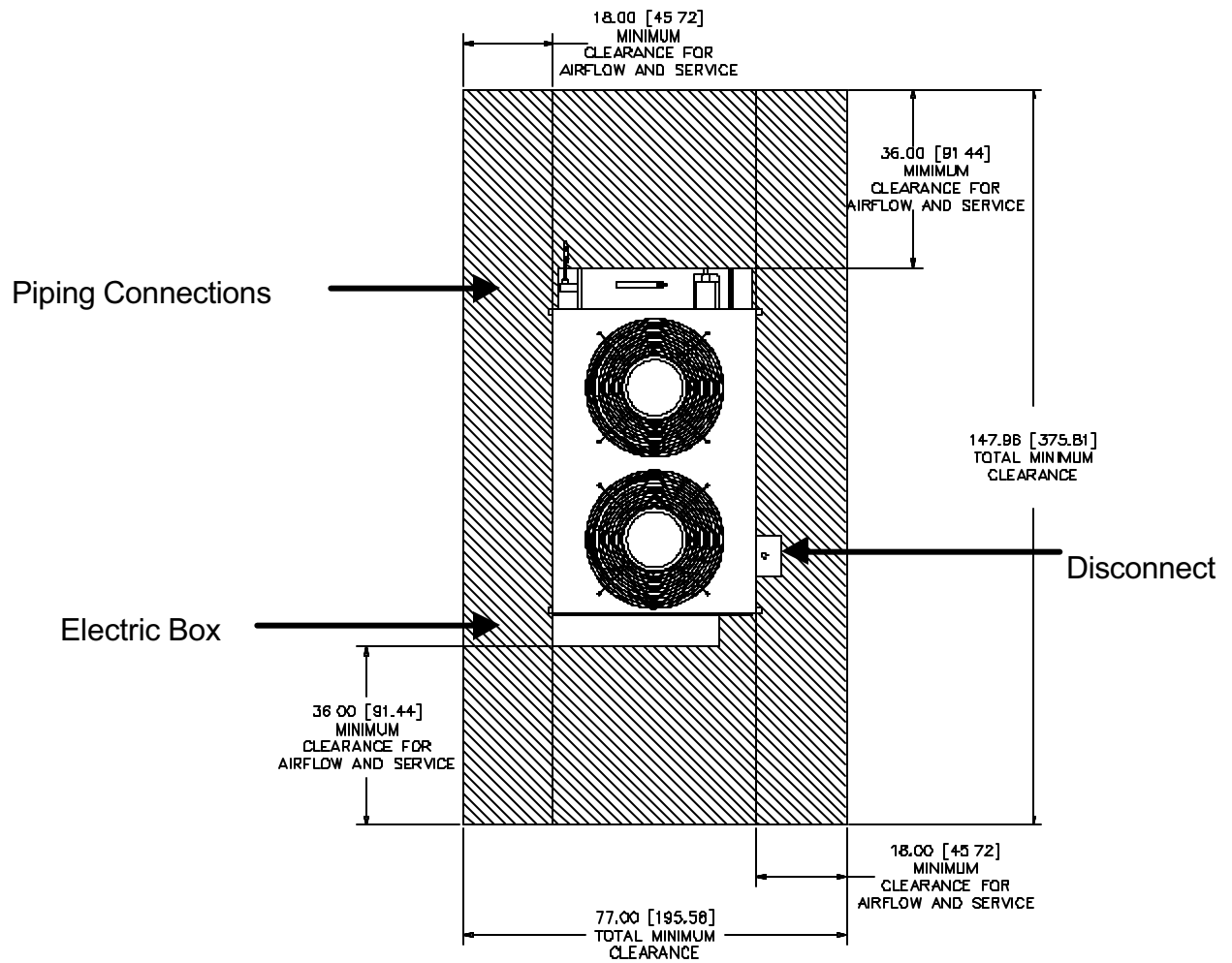
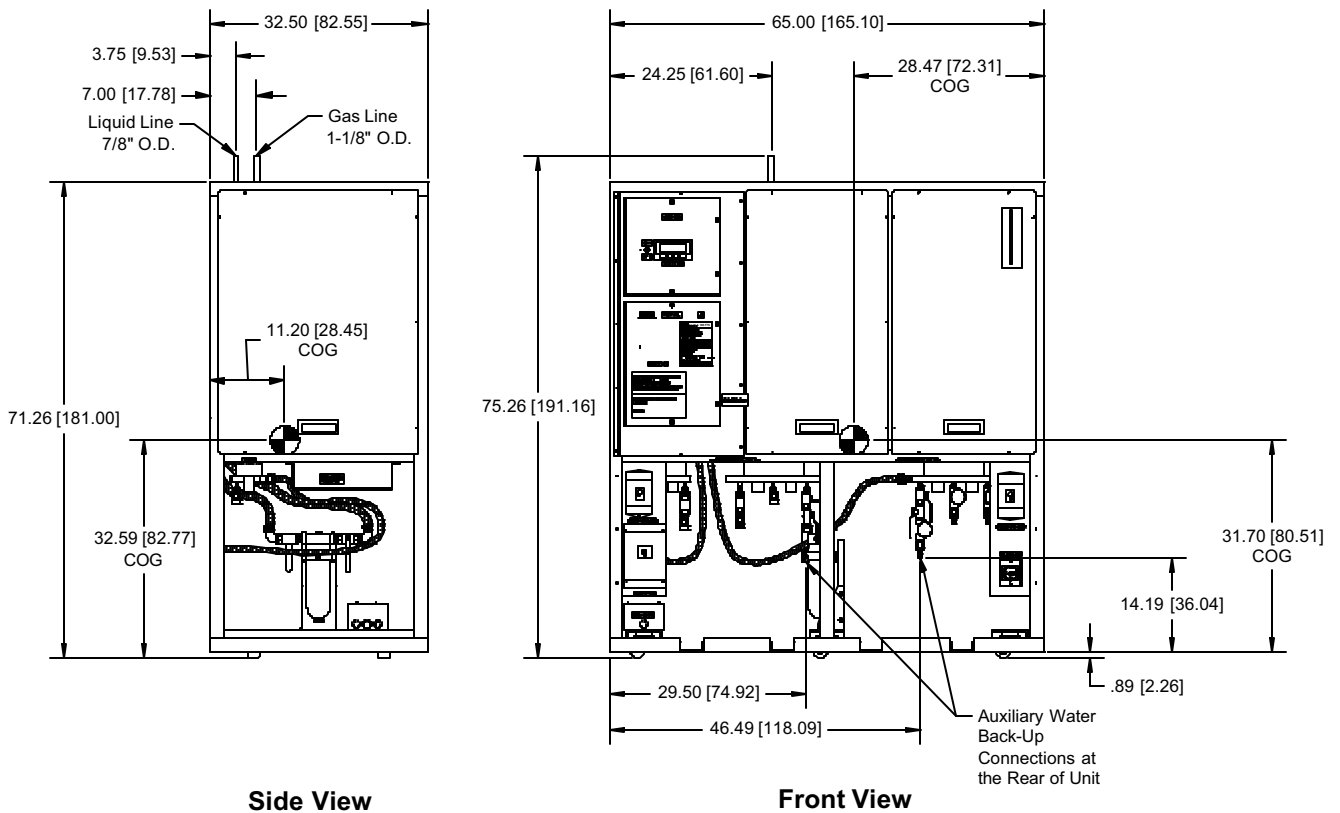


Illustration 5 **NECESSARY CLEARANCES FOR RAC ([] = cm)**

NOTE: RAC cannot be mounted with the airflow in a horizontal direction. Service area must be kept clear for proper airflow.

3-4 COMPONENT DIMENSIONS

Component	Height	Width	Depth	Installed Weight
SCC1	75.26 inches (191.16 cm)	65 inches (165.10 cm)	32.5 inches (82.55 cm)	1770 pounds Note
SCC2	42.12 inches (107 cm)	34.12 inches (86.67 cm)	40.5 inches (102.87 cm)	550 pounds
RAC	44 inches (111.76 cm)	81.50 inches (207.01 cm)	48 inches (121.92 cm) Including disconnect	470 pounds

NOTE: With Shield Cooler and Cryo Cooler installed.Table 7 **Uncrated Dimensions**Illustration 6 **SCC1 OVERALL DIMENSIONS ([] = cm)**

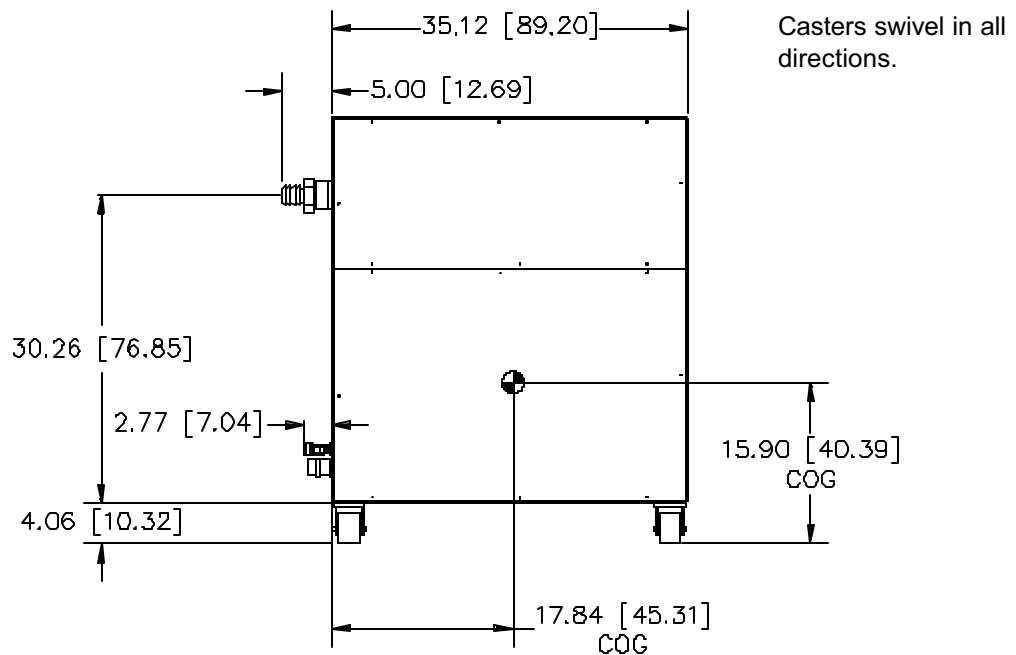
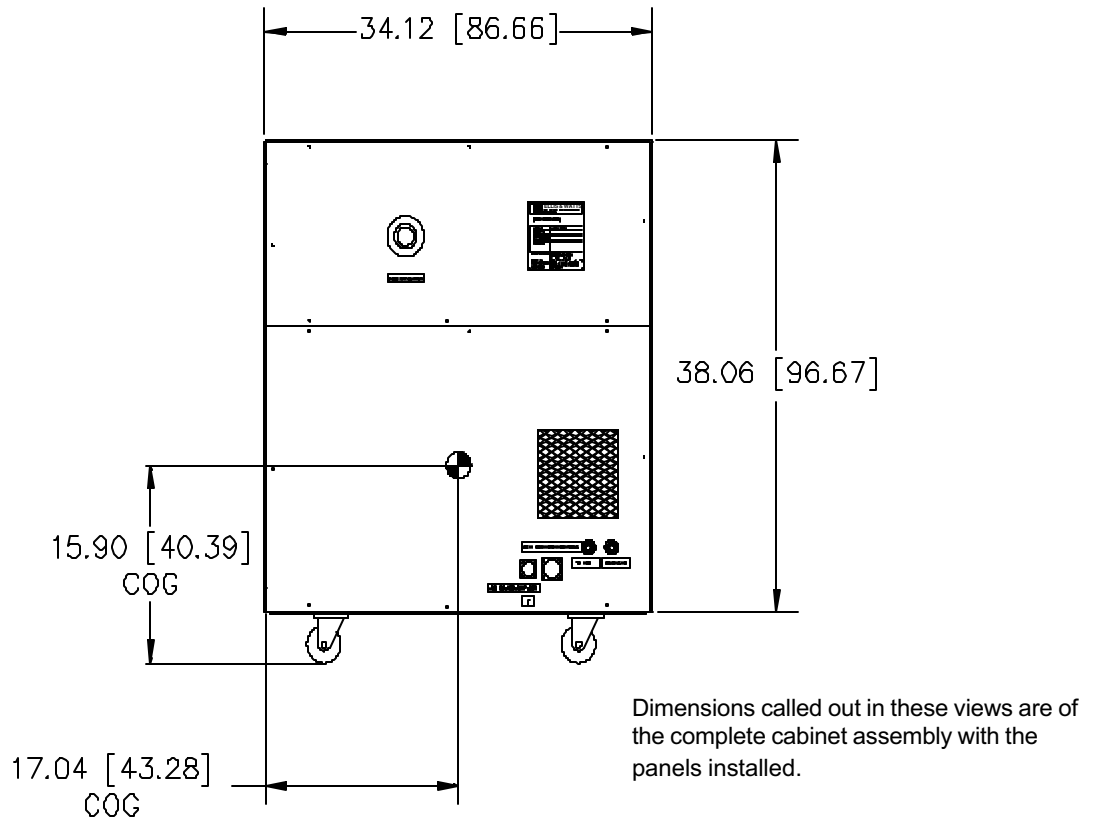


Illustration 7 **SCC2 OVERALL DIMENSIONS ([] = cm)**

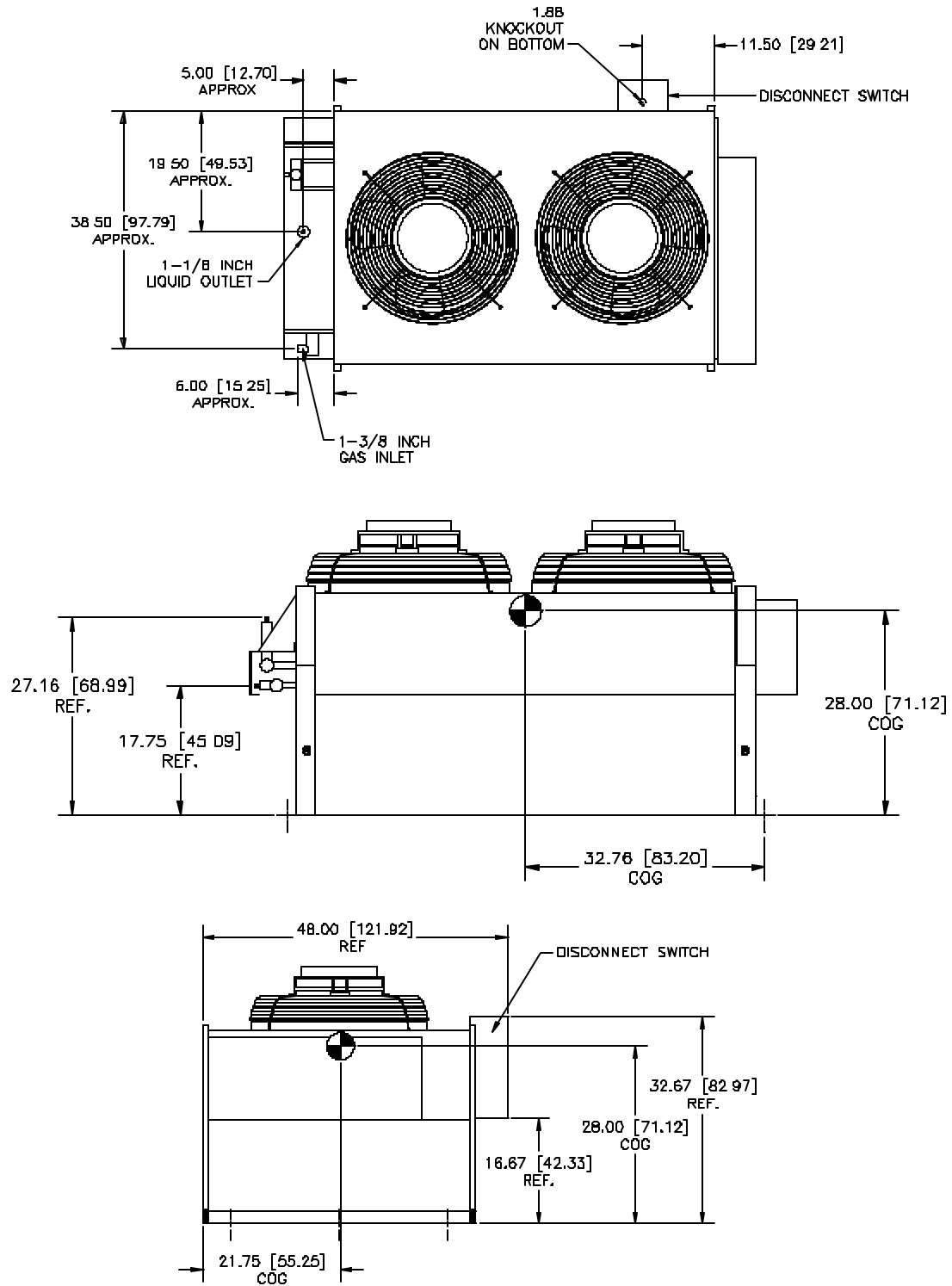
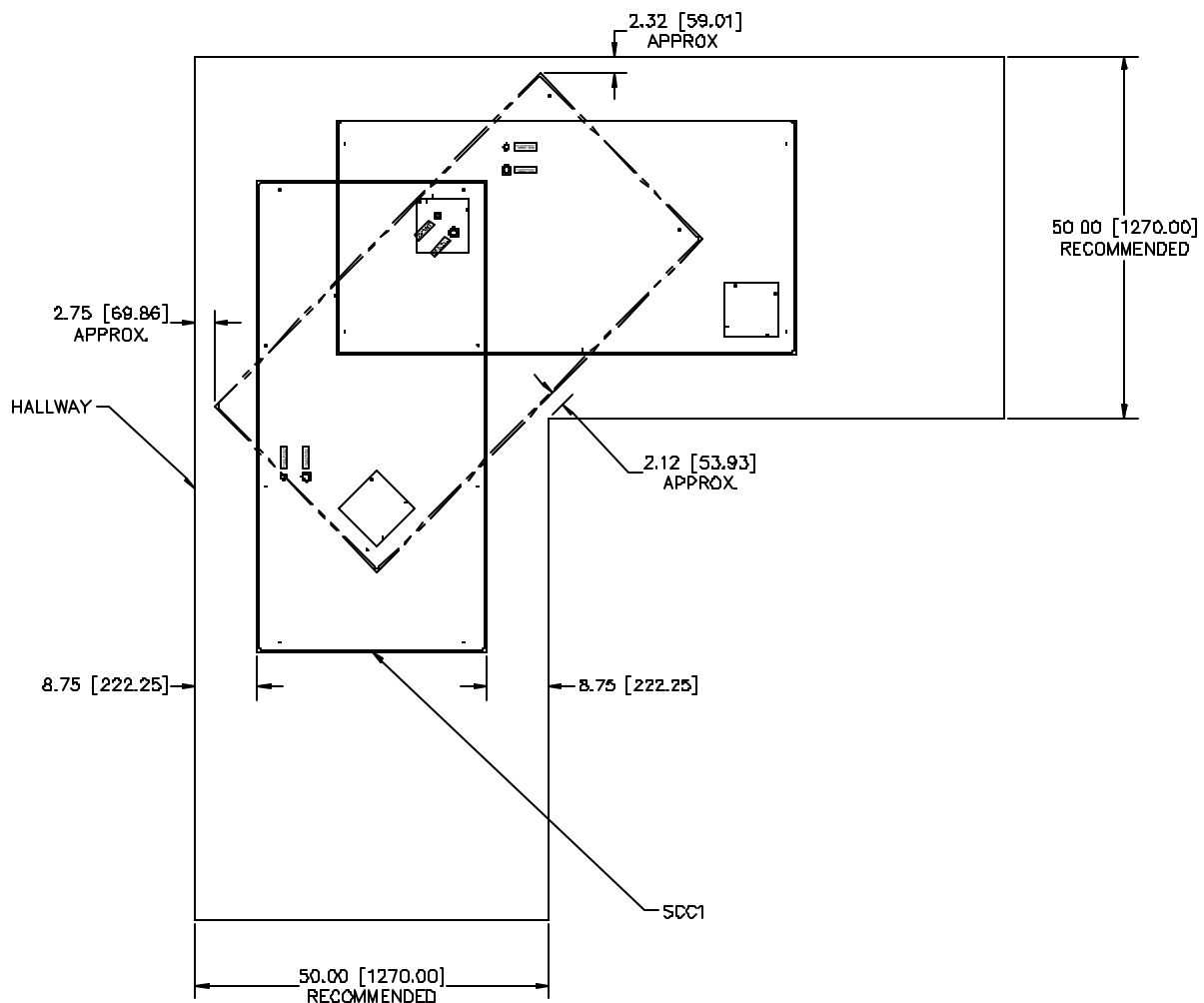


Illustration 8 **RAC OVERALL DIMENSIONS ([] = cm)**

3-5 MINIMUM DOOR/HALLWAY SIZES

Component	Hallway / Door Height	Hallway Width NOTE	Standard Doorway Width
SCC1	76 inches 193.04 cm	50 inches 127 cm	40 inches 101.60 cm
SCC2	43 inches 109.22 cm	42 inches 106.68 cm	36 inches 91.44
NOTE: Minimum width for 90° turns			

Table 8 **Door / Hallway Sizes**Illustration 9 **MINIMUM HALLWAY DIMENSION FOR TURNS ([] = mm)**

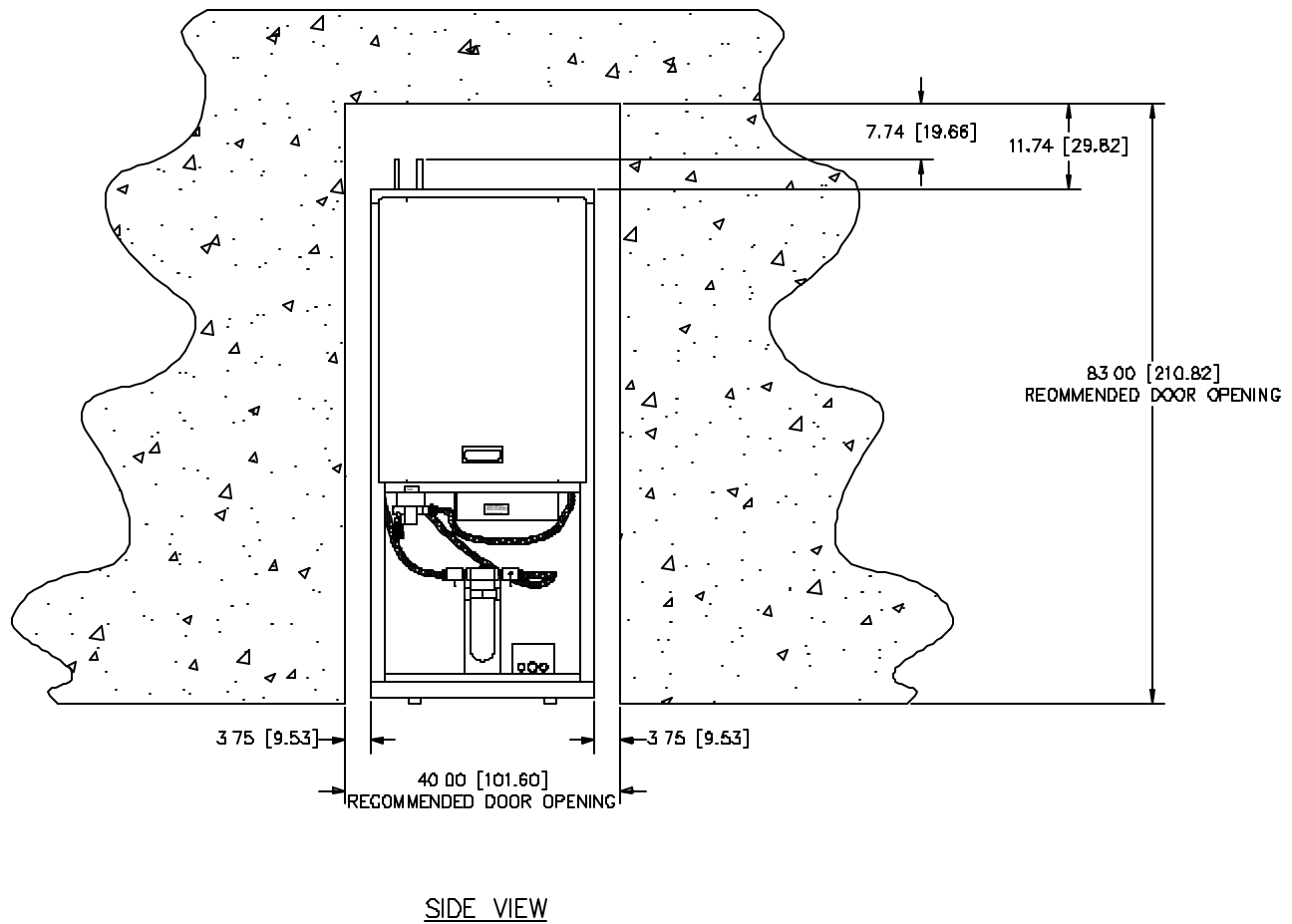


Illustration 10 **MINIMUM DOORWAY DIMENSIONS ([] = cm)**

3-6 DOORWAY TREASHOLD CONSIDERATIONS

The clearance from the bottom of the SCC1 cabinet to the floor is approximately 0.89 inches (2.26 cm) therefore doorway threshold should be as smooth as possible or special considerations have to be made to lift cabinet over a threshold height of 0.5 inches (1.27 cm).