

SECTION 5

MOVING MAGNET TO MR SUITE

Note

Roller dollies are recommended for moving magnet into the magnet room as shown in Illustration 5–1. If roller dollies are used, place steel floor plates along the magnet delivery route.



Use the jacking angles located on the magnet to raise the magnet for roller dolly installation.

Centerline of jacks must be at least 2" (50.8 mm) from end of angle and centered on horizontal leg of angle.

Use shims when rolling magnet with roller dolly over door thresholds and other inclines.

Push and pull magnet by the feet only.

Note

The magnet location in the magnet room must be marked. Use tape to mark the four corners of magnet on the magnet room floor. Refer to architectural drawings to determine the exact location of magnet within the magnet room.

1. Magnet is shipped with jacking angles installed. See Illustrations 5–1 and 5–2.
2. Push magnet to magnet room. If using a motorized tow vehicle, attach chains only around the magnet feet. Apply no loads to any other parts.
3. If there are turns in the delivery route, adjust roller dollies to appropriate positions to negotiate turns.
4. Move magnet into the magnet room. If there is not enough vertical clearance for the magnet, call the MAC Team representative for removal of vertical penetration components.

Note


Table end of magnet must be oriented away from rear wall of exam room.



Keep magnet level at all times. Uneven jacking of corners could result in shifting of magnet on jacks.

5. Jack the magnet up sufficiently, at the 4 jacking brackets, and remove the 4 roller dollies.
6. Slowly lower magnet onto the magnet room floor. Release pressure simultaneously in pair of jacks on one end until 1 to 2 inches (25.4 mm to 50.8 mm) lower than opposite end. Repeat end to end until feet are on floor.

ALL DIMENSIONS ARE IN INCHES (MILLIMETERS)

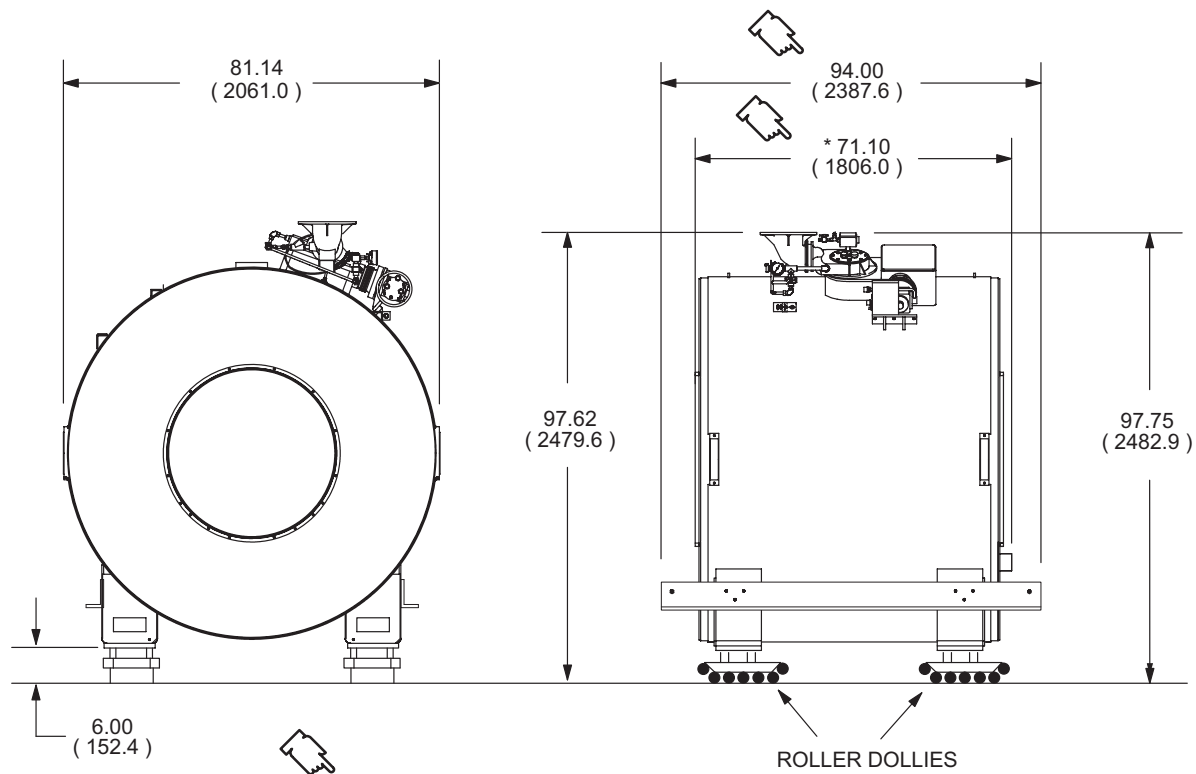

* INCLUDES BRM MOUNTING BRACKET AND HEX NUTS


MAGNET WITH ROLLER DOLLIES ATTACHED

ILLUSTRATION 5-1

Note

The height can be reduced by putting roller dollies under the lifting rails.

ALL DIMENSIONS ARE IN INCHES (MILLIMETERS)

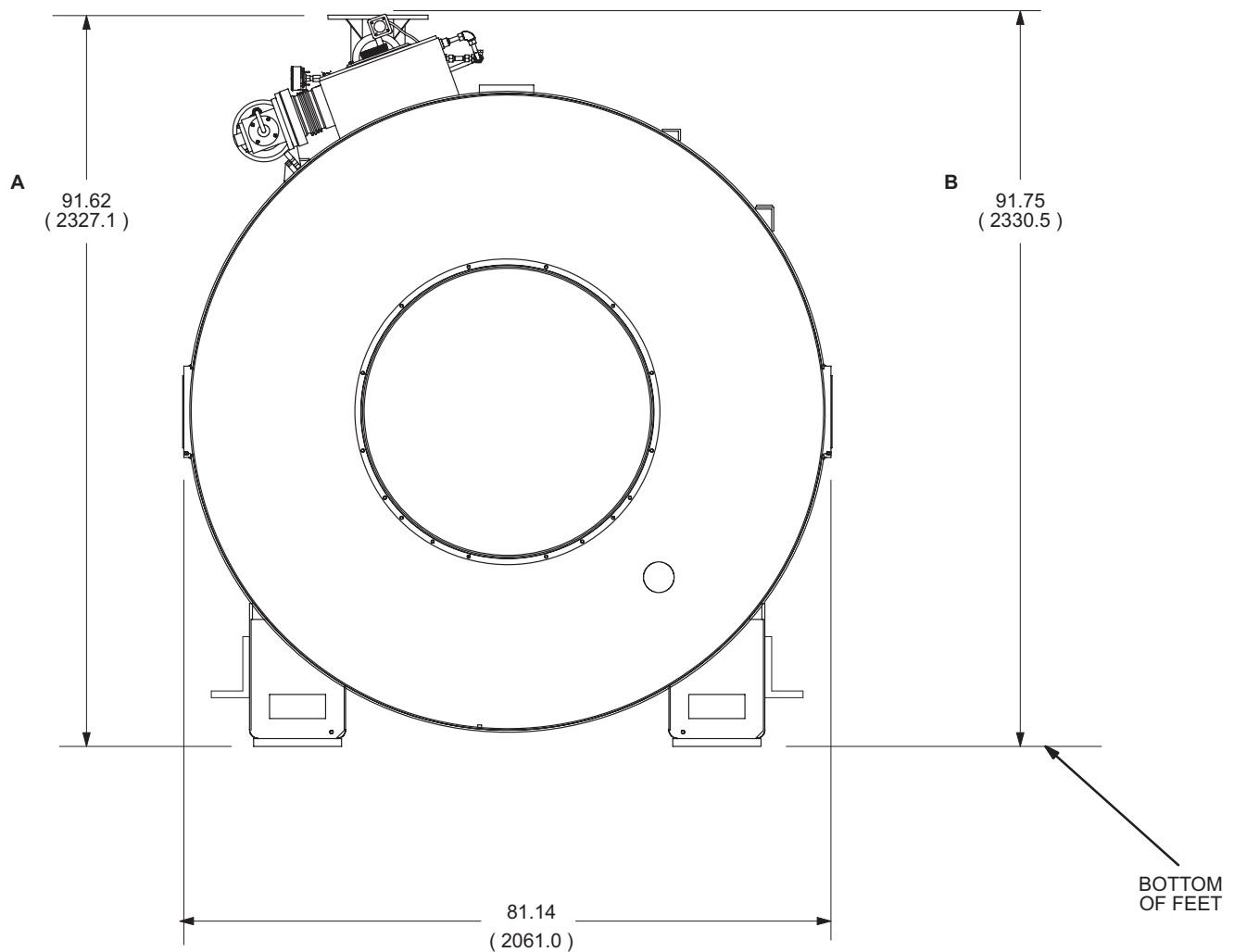

A: DISTANCE TO TOP OF VENT ADAPTER
B: DISTANCE TO TOP OF SHIM LEAD CONNECTOR BOX

Note

Minimum height for service clearance is 101.86" (2587.2 mm). This clearance is needed for ramp lead, shim lead, and fill line installation.

MAGNET CLEARANCE DIMENSIONS
ILLUSTRATION 5-2