

MR Monthly CSE-TSE Conference Call, September 2018

Jerry Howze, Director Service Support



Restricted © Siemens Medical Solutions USA, Inc., 2016

Page 1

MR Monthly CSE/TSE Conference Calls

Objectives



Agenda:

- Safety First!
- MR – Quality Initiatives
- Updates
- Speed Infos
- Tips of the month
- Max Update
- Special Topic
- Feedback from field

- Objective of monthly conference calls
 - Build stronger connection between the Cary campus and the field
 - Provide useful/timely information to our field CSEs
 - Share information on improvement projects
 - Gather feedback from the field

EHS Topic Brief Ladder Safety

Introduction: Joe Fass



Joined Siemens October 2015 – Cary, NC

- Service organization EHS

Background

- Operations – Mfg./Industrial Engineering
- Environmental, Health and Safety

My Approach

- Make something better today than it was yesterday
 - $\frac{3}{4}$ Don't let perfect get in the way of better
- Safety is what we do, not what we avoided

Personal

- Raleigh, NC with my wife and kids
- Love the outdoors and outdoor activities, especially fishing



Ladder Safety – 10 tips to keep you safe

1. Carefully inspect the ladder for defects, checking for cracks, corrosion and that bolts and rivets are secure. Tag and remove unsafe ladders from service.
2. Make sure the ladder's feet work properly and have slip-resistant pads.
3. Use a fiberglass ladder (electrical safety)
 - MRI Platform (secured) is preferred for most tasks
4. When setting the ladder, look for a safe location with firm, level footing and rigid support for the top of the ladder. Be sure to set it at an angle per the manufacturer's guidance (fully extend).
5. When climbing off a ladder at an upper level, make sure the ladder extends 3 feet above the landing (extension ladders).



Ladder Safety

6. Use three points of contact — keep 1 hand and both feet or both hands and 1 foot in contact with the ladder at all times.
 - Climbing: 2 feet, 1 hand or 2 hands, 1 foot
 - NOTE: When necessary for two handed task: 2 feet and waist/torso leaning against ladder is acceptable
7. Never carry any load that could cause you to lose balance.
8. Never stand on top of a ladder.
9. Don't pull, lean, stretch or make sudden moves on a ladder that could cause it to tip over.
 - $\frac{3}{4}$ Basically keep your belt buckle between the rungs
10. Avoid setting the ladder near exit doors, near the path of pedestrian or vehicular traffic.



Environmental Health and Safety Discussion/Questions - Contact Info

Joe Fass
Environmental, Health and Safety
HC NAM USA - Service

221 Gregson Drive
Cary, NC 27511, USA

P: +1 919-461-6020

M: +1 919-239-5102

Joseph.Fass@siemens.com



Agenda:

- MR – Quality Initiatives

MR Quality & Productivity Initiatives

Kyle Hamilton

Kyle.hamilton@siemens-healthineers.com


Strategic Business Planner
Strategy and Projects

MR Quality Initiatives

Material Consumption (Contract MS, MM notifications only)

Agenda:

- MR – Quality Initiatives

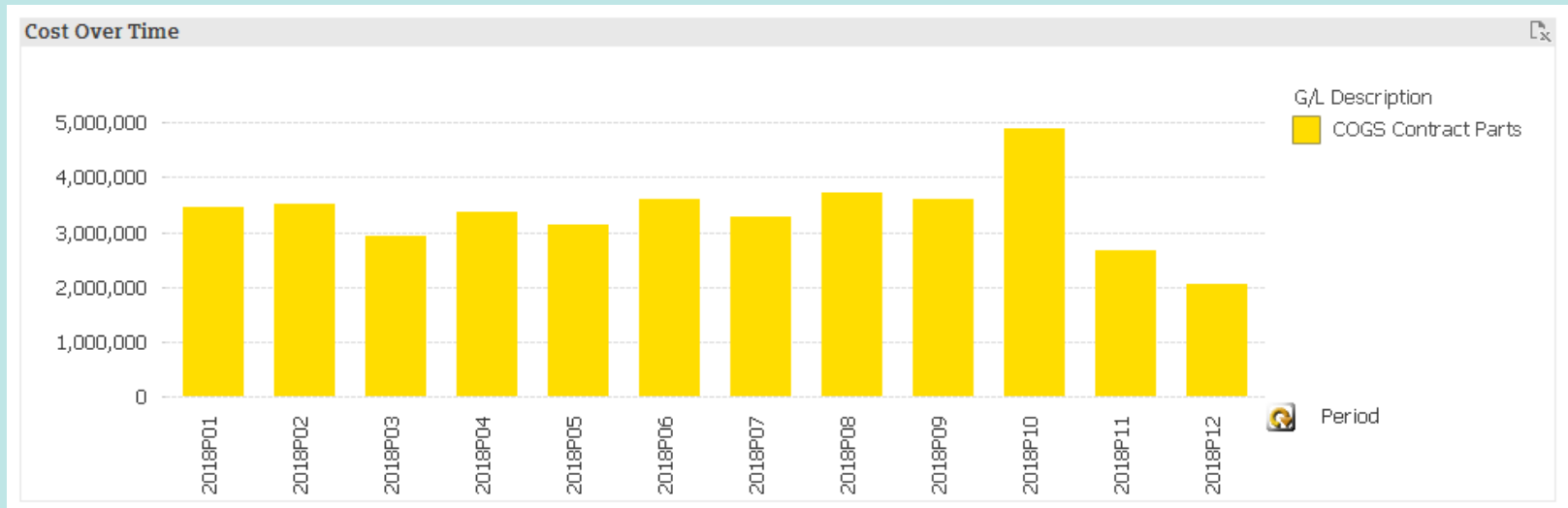


MR Quality Initiatives

Field Labor (Contract MS, MM notifications only and regular hours)

Agenda:

- MR – Quality Initiatives

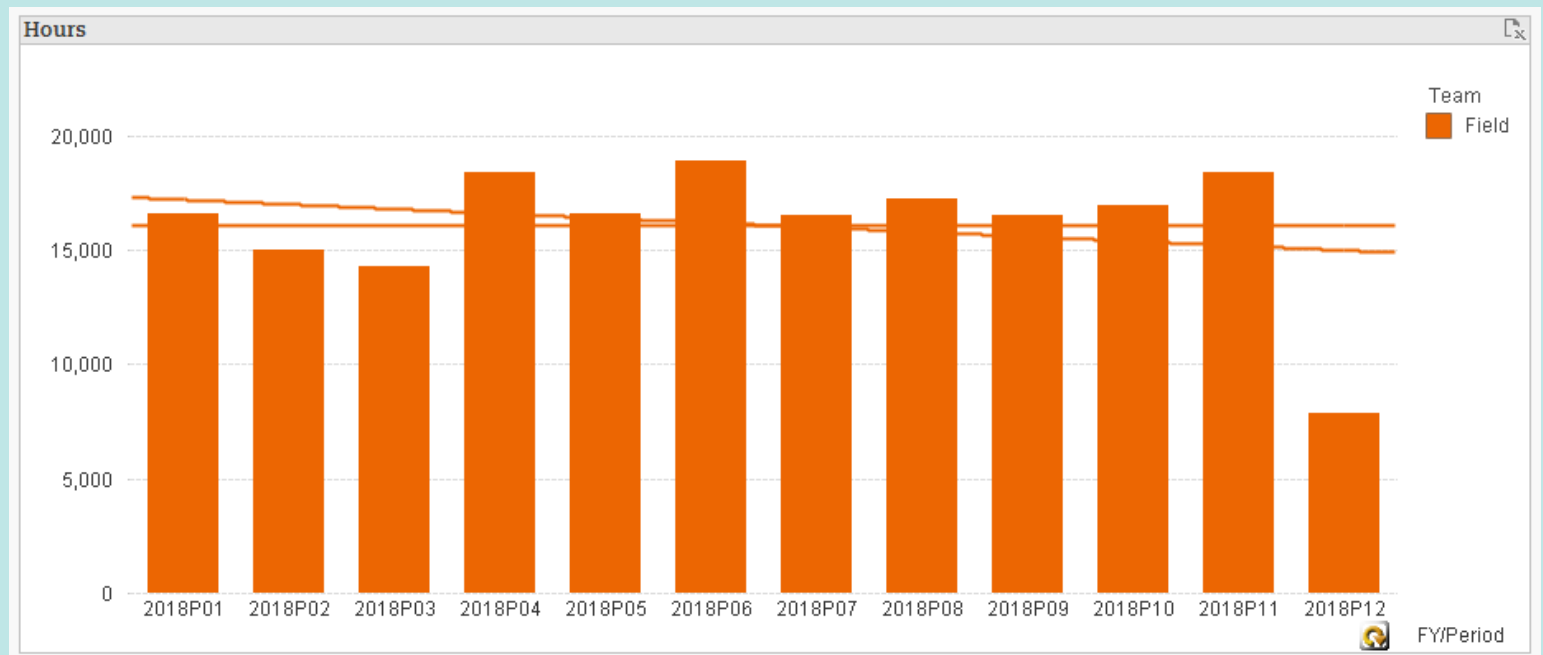


MR Quality Initiatives

Helium Spend

Agenda:

- MR – Quality Initiatives

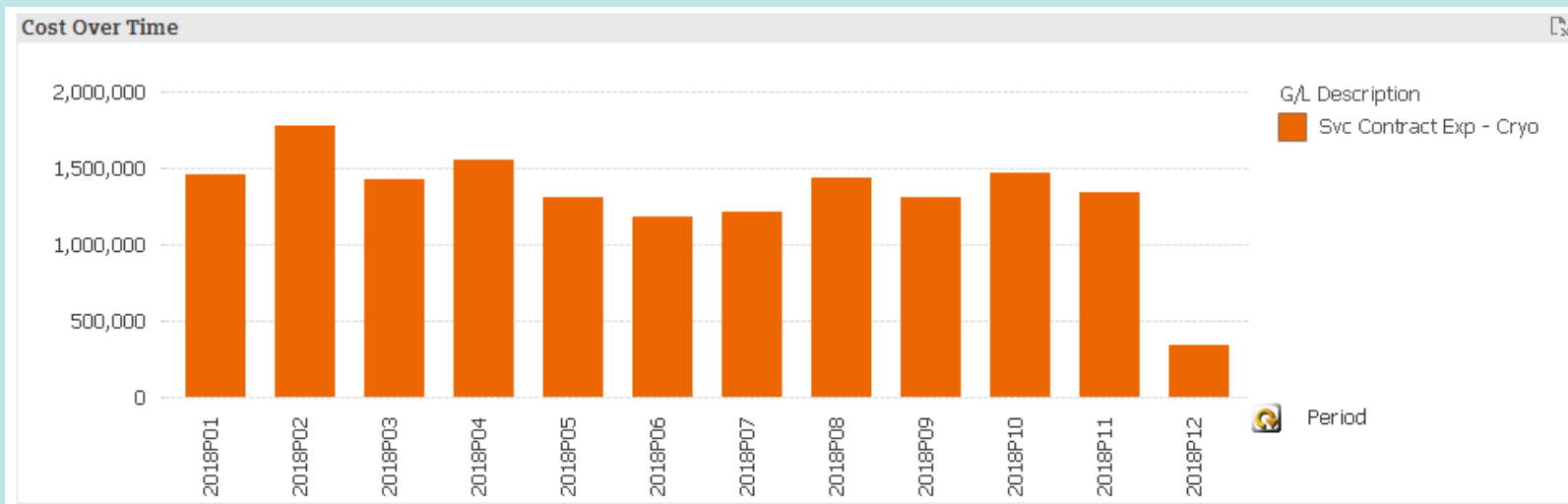

	FY17 Liters	FY18 Liters	Reduction
Q1	568,227	447,520	120,707
Q2	598,922	582,230	16,692
Q3	547,152	569,530	-22,405
YTD	2,050,915	1,984,442	66,473


MR Quality Initiatives

MR Roadshow

Agenda:

- MR – Quality Initiatives

- Cary, NC
- NYC
- Minneapolis, MN
- Boston, MA
- Atlanta, GA
- Phoenix, AZ
- Tampa, FL
- Los Angeles, CA
- Dallas, TX
- San Francisco, CA
- Houston, TX
- Seattle, WA
- Columbus, OH
- Baltimore, MD
- Chicago, IL
- St. Louis, MO

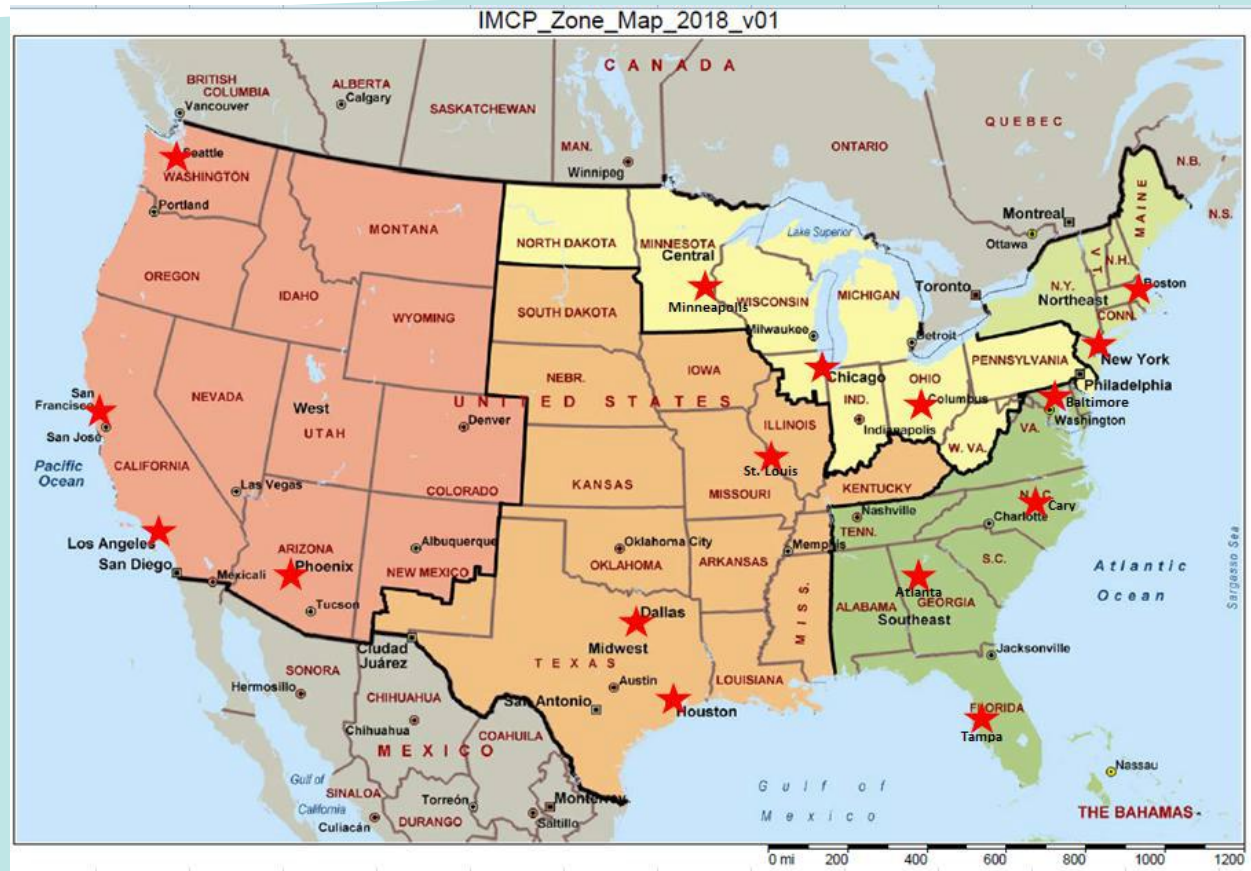

MR Quality Initiatives

MR Roadshow



Agenda:

- MR – Quality Initiatives

	Dates	Location
Week 1	Oct 4th and 5th	Cary, NC
Week 2	Oct 8th and 9th	NYC
	Oct 10th and 11th	Minneapolis, MN
Week 3	Oct 15th and 16th	Boston, MA
		Atlanta, GA
	Oct 17th and 18th	Phoenix, AZ
		Tampa, FL
Week 4	Oct 22nd and 23rd	Los Angeles, CA
		Dallas, TX
	Oct 24th and 25th	San Francisco, CA
		Houston, TX
Week 5	Oct 29th and 30th	Seattle, WA
		Columbus, OH
	Oct 31st and Nov 1st	Baltimore, MD
		Chicago, IL
		St. Louis, MO

Topics

- Local Coils
- Guardian
- Vida
- Spike Troubleshooting
- Tables
- General Troubleshooting
- Onsite to Online
- Miscellaneous Topics
- Open Q&A

MR Quality Initiatives

IQ Team and Coil Reviews



Agenda:

- MR – Quality Initiatives

- Automated emails when specific coils are ordered
 - Spine Matrix, Head Matrix, Body Matrix, Body 18
- Research notification and flag based on 3 scenarios
 - Coil replaced defective
 - Coil replaced investigation
 - Coil replaced IQ
- 15 out of 33 coils ordered since 9/6 marked for monitoring
 - 7 flagged for IQ review
- Taking a more proactive approach to ongoing coil issues at a site

Agenda:

- Safety First!
- MR – Quality Initiatives
- Updates
- Speed Infos
- Tips of the month
- Max Update
- Special Topic
- Feedback from field

MR Technical Updates

Updates, SKBs (SI), Tips of the month

MR Technical Updates

Updates, Speed Infos, Tips



Agenda:

- Safety First!
- MR – Quality Initiatives
- **Updates**
- Speed Infos
- Tips of the month
- Max Update
- Feedback from field

- *MR006/18/P-Replacement of string encoder. The string encoder could fail and the tabletop could not be moved by motor. (Vida)*
- *MR008/18/P-Replacement of CFDA label. A new label to replace old CFDA label which was printed and placed on the delivery systems. (Essenza Dot)*
- *MR011/18/P-Mismatch between MR and PET images. This Update Instruction is to inform our mMR customer that we identified a mismatch of PET and MR images in z-direction when viewing fused images on variant mMR system. (Biograph mMR)*
- *MR039/17/P-Magnetic Handles of CDT. The springs, screws and threaded inserts in both handles of the transfer boards are magnetic and can cause signal loss in the region of knee and foot. (Mobile Table-Aera, Skyra)*
- *MR017/18/R-- Fail of Auto-Restart Chiller Pump. To improve the reliability of the automatic restart of KKT Chillers after a power failure, the pump configuration will be reconfigured. (Aera, Skyra, Verio, Avanto, Espree, Av. Fit, Verio Dot, Prisma, Prisma Fit, Skyra Fit, Vida)*
- *MR021/18/P-- Removal of tank anode rod in the chiller. (Amira, Spectra)*
- *MR009/18/R-Wrong label on 75/22 Power Stages (Skyra)*
- *MR013/18/P-N4 VD14A-SH07 Software update (Essenza Dot)*
- *MR019/18/P-VE11S-SP02 Software update (Amira, Sempira)*

MR Technical Updates

Updates, Speed Infos, Tips



Agenda:

- Safety First!
- **MR – Quality Initiatives**
- Updates
- **Speed Infos**
- Tips of the month
- Max Update
- Feedback from field

- SKB0097288-Additional GC transport pallet is required for a GC exchange. The spare part Gradient Coil for Skyra systems (part-no. 10358500) or Gradient Coil for Vida systems (part-no. 10678889) is now delivered on an improved GC exchange palette. (Skyra, Vida)
- SKB0093764-Low spikes on Z axis. (Aera XQ, Vida)
- SKB0090932-No refund for returned RF amplifiers TBX or TX-Boxes without filled in questionnaires. To get a refund, return the replaced RF amplifier TBX or TX-Box with a completely filled in questionnaire. Link to questionnaire:
 - <https://bps-healthcare.siemens.com/productinformation/ProductInformation.aspx?uuid=7abf6bb6-5daf-440c-9920-fab2d4640084>
- SKB0093336-Online Help for TU step BC Coil Power Losses is wrong. VD14A, VE11Q (Essenza Dot)
- SKB0097973-MR View & Go - SW User Interface- SW crash- "Blue-screen" after pulling out USB stick. Unintended switch on/off by user.(Vida)
- SKB0089091-Checking the EPO circuit. (MR015/17/P) (MR)
- SKB0098016-Abdomen DOT Engine TB- wrong timing in the dynamics. (Vida)
- SKB0084468--Buzzing sound of Operation Unit (p.n. 11060844) due to a new version of the Operating Unit Power Supply. (Aera, Skyra, Aera24, Av. Fit, Prisma, Prisma Fit, Skyra Fit, Spectra, Terra)
- SKB0075849-New spare part for Patient Table and Table tops with wide notches. Two types of notches (6mm & 10.5mm)
- SKB0098620-PET-service within SESO does not start after SD05 install. (Biograph mMR)
- Check with DVM and adjust voltages. (Vida)

MR Technical Updates

Updates, Speed Infos, Tips

Agenda:

- Safety First!
- MR – Quality Initiatives
- Updates
- Speed Infos
- Tips of the month
- Max Update
- Feedback from field

- SKB0040078-Handling of Li-Battery Pack (#10128924). The new Li battery pack is used with MSUP 10523587 and 10892850. It is meant to provide an energy backup for quenching the magnet in case mains power fails. (MR)
- SKB0043618-Ramp Currents. Ramping the magnet without the acceptance documentation. (MR)
- SKB0060368-Intermit. Communication Error with Patient Table. On Avanto Dot and Verio Dot systems intermittent communication problem between host and patient table has been reported.
- SKB0085412-Local and Global DIXON Fat-Water Swaps. Some cases it could happen that the system exchanged the fat and the water signal. (MR)
- SKB0096411-Reducing Diffusion Artifacts with DWI Resolve & MDDW sequences.(Amira, Spectra)
- SKB0099295-Remote updates will NOT work in case VB19A-SH08.(Verio,SATS,TATS,Avanto, Esp)
- SKB0099082-- When is it necessary to change the ERDU battery? The ERDU backup battery 10128925 (#10128924) is built from Li metal batteries, is not rechargeable and has an estimated life time of more than 10 years. This is valid as long as it is not getting unloaded or discharged, i.e. a quench button is pressed in a mains power failure scenario. (MR)
- SKB0098949-Exchange of MSUP (#10523587 or 10892850) after a quench. We recommend **not** to swap the MSUP, the alarm box, the ERDU battery or the DAFID after a quench. (MR)
- SKB0098940-- Proper installation of the clamps for the gradient cables. Applying too much torque can break the clamps for the gradient cables and can also damage the gradient cables themselves. (MR)
- SKB0085720-- Reconstruction error. With spectroscopy when multiple 'Measurements' or 'Save single averages' are enabled together with 'Prescan Normalize'. (Aera,Aera24,Av. Fit,Prisma,Prisma Fit, Skyra, Skyra Fit)
- SKB0100133-- Seso cannot be opened as syadmin after installation of MR025/17/R (Aera,Skyra,Aera24,Av. Fit,Prisma, Prisma Fit,Skyra Fit)
- SKB0100610-- New Tune up step "Frequency Calibration" with VA11A. (Sola,Vida)

MR Technical Updates

Updates, Speed Infos, Tips

Agenda:

- Safety First!
- MR – Quality Initiatives
- Updates
- Speed Infos
- **Tips of the month**
- Max Update
- Feedback from field

- Introducing the Known Issue.
- A known issue is the repeated observation of a possible defect, involving a system, software, or reagent, which:
 1. Has occurred on 3 or more systems
 2. Has **NO** root cause determined by PMS (Problem Management)
 3. Is without an SKB support document already
 4. Does not have a workaround / interim solution
- Why is it important?
 - A known issue is important because it bridges the communication gap from the time an issue is discovered
 - until a workaround or final solution is communicated in the form of an SKB. A known issue may be seen
 - during new product introduction or after new UI's have been released. Knowing about a Known Issue can
 - help save hours and parts for issues that have no solution and are still under investigation.
- Where is the value?
 - By being aware of a known issue, you can reduce the amount of site visits, parts usage, and increase
 - customer transparency, ensuring that we're getting the right solution to the system as opposed to just
 - another part.

MR Technical Updates

Updates, Speed Infos, Tips

Agenda:

- Safety First!
- MR – Quality Initiatives
- Updates
- Speed Infos
- **Tips of the month**
- Max Update
- Feedback from field

- Ticket Escalation
 - If you need support because your customer has an issue with our Siemens system, we want to help this customer so that this issue is solved as soon as possible and we don't want to waste time with un-necessary actions. Depending on the issue we need specific data and information to understand the issue. Save-logs, DICOM images, Camtasia video, Screen-shots and log-files are just a few ways of getting information. Please see SKB0097975 for further details.
1. System Type? SW version? Serial Number?
 2. When did the issue happen the first time? (Date and Time)
 3. What exactly is the problem?
 4. Is there an error message?
 5. What happened beforehand? (SW-Update, SW crash, protocol change)
 6. Can you reproduce the issue or is it intermittent?
 7. Is it getting worse?
 8. How often does it happen?
 9. Was Siemens protocol/workflow/DOT used?
 10. Any other system problems that could be related? (Coils, Adj. issues)

Chiller Support

Eugene Murano/Tridib Borkakoty

Agenda:

- MR – Chiller Support
- MR – Quality Initiatives
- Updates
- Speed Infos
- Tips of the month
- ERDU Update
- Quench Vent
- PM Transition
- Max Update
- Feedback from field

- **Water Leaks**
- **UIs for Chiller**
- **Tips**
- **MR Chiller Support Team**

- **Water Leak – Who is responsible?**
 - Within the Chiller - KKT/GBT
 - Piping between Chiller and IFP (including connections at both ends) - Customer
 - IFP and within our equipment - Siemens
 - Top-offs by KKT/GBT are billable if the glycol loss resulted from a confirmed leak outside the chiller.
- **Glycol Top-Off (ECO Chiller)– Should CSE perform top-off?**
 - If glycol is lost when replacing a water cooled component and auto-fill device didn't re-pressurized the system
 - CSE should avoid topping off, especially ECO chillers, when the reason for glycol loss is unknown. In such cases CSE should perform a general check of the system for any leak in MR system and escalate it to tech-support if KKT/GBT techs are also unable to find any leaks. We will work with KKT/GBT to determine next steps.
 - ECO chillers must be turned-off during top-off. The static pressure should be read on the IFP gauge and adjusted according to CBDOC values.

- **UI MR034/17/P for T Fitting- Cancelled**
 - UI cancelled because of defective component (Elbow) in the kit
 - T-Fittings and Elbow cracks will be now repaired reactively (only when it fails) free of cost by KKT.
 - We are reducing risk of T-Fitting/Elbow cracks by configuring ECO133 chiller as ECO122 for sites having XJ gradients (45kW)
 - ECO133 can be configured as ECO122 by just pulling out a jumper
 - The small compressor will be disabled when configured as ECO122
- **MR017/18/R FAIL OF AUTO-RESTART CHILLER PUMP**
 - ECO Chillers Manufactured between June 2014 and Sept 2017.
 - The wrong Pump setting causes a 45 minutes delay to auto reset the Pump VFD after a phase imbalance or loss.
 - The UI will be performed by SBT/KKT at sites having chiller contract during next PM.

- **KSC215 Night Mode (FW Version 3.1 – 3.3)**

- In night mode, when MR system is in Standby or OFF, the chiller water temperature allowed to vary from 20 to 26 degree Celsius
- Unlike ECO Night mode, the pump speed in KSC215 is fixed

- **ECO Chiller Fan Error**

- If the chiller shuts off intermittently and/or MR system generates a FAN1 Fault (Frequency Converter error), it is very likely that the chiller has a defective fan and it will be completely down if not replaced soon.
- IFP LED Status

LED Line PCB-Board	1	2	3	4	5	6	7	8
Active LED's			X	X	X	X	X	X

- **Defective IFP Pressure Sensors**

- MR CSE is responsible for IFP component replacements.
- No need to drain glycol while replacing the IFP pressure sensors.
- Sometimes IFP sensors read abnormal values if the bottom of the sensor inside the piping (the stator core) get clogged with dirt or debris. The core can be removed, cleaned, and installed back using a tool available in refrigeration stores without draining the IFP. The Sensor replacement process is attached.

IFP TYPES

- IFP 098



- IFP 084



- IFP 090



Fault you will see in Event log thru SRS

This Fault Occurs when there is a blockage underneath sensor or sensor is defective showing a low pressure reading.

Remedy that has worked so far :

- 1.Replacing stator core under transducer
2. Replacing sensors

LED Line PCB-Board	1	2	3	4	5	6	7	8
Active LED's		x						x

Pipes Connected Wrong or Sensor on IFP Changed

LED 2&8 are activated, the difference between sensor 8B2 IFP and 8B3 IFP
is negative

See also manual page 60


Tool Needed for Stator core Removal.

Both can be found at any AC or Refrigeration supply house

Option #1 allows you to remove stator core with out draining system




Option #2 will need to drain system with this tool

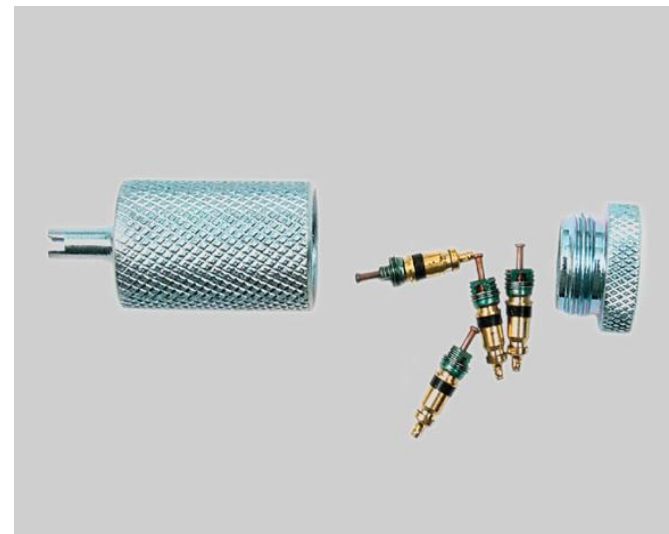

Location of Pressure Transducers


- 1- Supply Pressure
- 2- Return Pressure

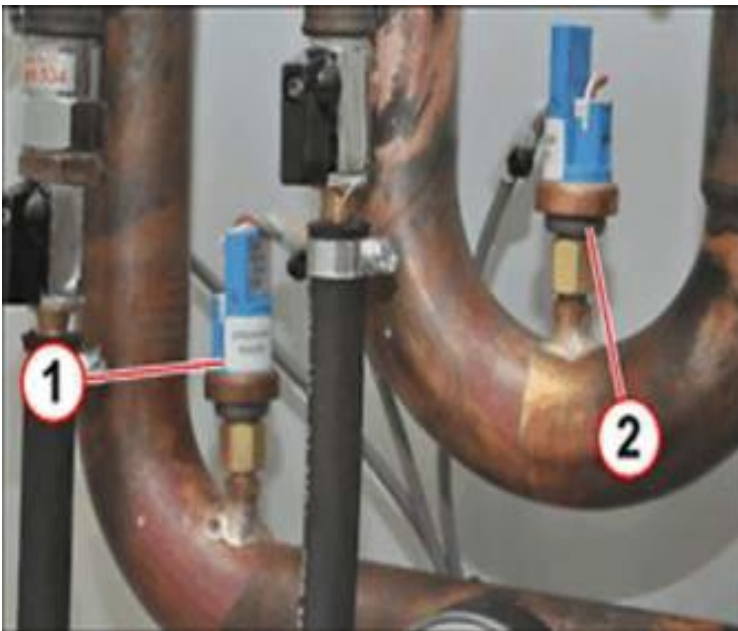

How to replace Also in CB DOCS

4.5.2 Preparatory Work Steps

- Switch the chiller to service mode
 - Switch off the compressor.
 - Switch off the chiller unit, 4Q1.
- Fig. 122: Chiller electronic cabinet

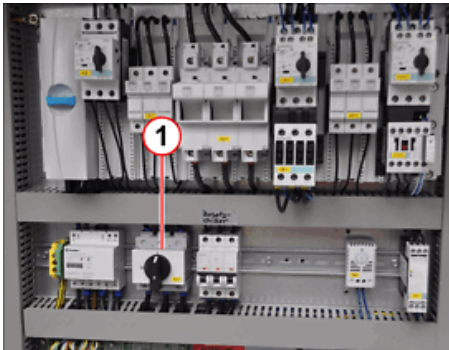

- Switch off EPC/F102 at the mains box
- Close the shutoff valves for the compressor

4.5.3 Removal

Water will not escape if the sensor is removed.

Loosen the screw of the sensor

Remove the sensor.

Disconnect the cable at terminal block X2
Supply water pressure sensor X2 - 16, 17

Return water pressure sensor X2 - 18,

4.5.4 Installation

- x Install the new pressure sensor.

Tighten the sensor with the screw

- x Connect the sensor cable to terminal block X2 - 16, 17 or at X2 - 18, 19

4.5.5 Startup

4.5.5.1 Mechanical Work Steps

1. Check the filling pressure
 1. Refill water if necessary
2. Switch on EPC/F102
3. Switch on the chiller 4Q1.
4. Switch on the compressor.

4.5.6 Final Steps

4.5.6.1 Final Checks

1. Check filling pressure after system is running
2. Refill glycol mixture (KKT Protect 37) if pressure dropped.
3. Check the filling pressure of the auto-refill device (if available).
4. Auto-refill device shall be filled at 6 bar
5. Refill glycol mixture (KKT Protect 37) if necessary.
6. Check the system for water leakage
7. Check the parameters of the cooling system via SeSo/Magnet & Cooling platform

Reset the chiller service mode

- **Functional Location Sticker for Chiller**

- Some chiller F/Ls are not properly associated with MR F/L and requires additional time to locate it.
- If customer can provide the chiller F/L, the chiller calls can be routed promptly and appropriately.
- Order Info for F/L stockers- Item# A91SVC-17941-P1-4A00
- <https://custompoint.rrd.com/xs2/prelogin?CMPID=3045&qwerty=18091707>

Agenda:

- **MR – Chiller Support**
- **MR – Quality Initiatives**
- **Updates**
- **Speed Infos**
- **Tips of the month**
- **ERDU Update**
- **Quench Vent**
- **PM Transition**
- **Max Update**
- **Feedback from field**

MR Chiller TSE

Eugene Murano

TSE MR Chiller Tech Support

221 Gregson Drive

27511 Cary, NC

919-319-2900 Ext 798804

919-883-8963 (mobile)

Questions?


SRS / CAN24

AJ Demianenko

- **Agenda:**
- Safety First!
- MR – Quality Initiatives
- Updates
- Speed Infos
- Tips of the month
- Max Update
- **Special Topic**
- Feedback from field

SRS/FTP/RUH/CAN24

- Lost momentum in recent months
- Encourage customers to reboot systems each week to maintain performance and establish a CAN24 Heartbeat check.
- CAN24 provides critical magnet data when the main system is off.

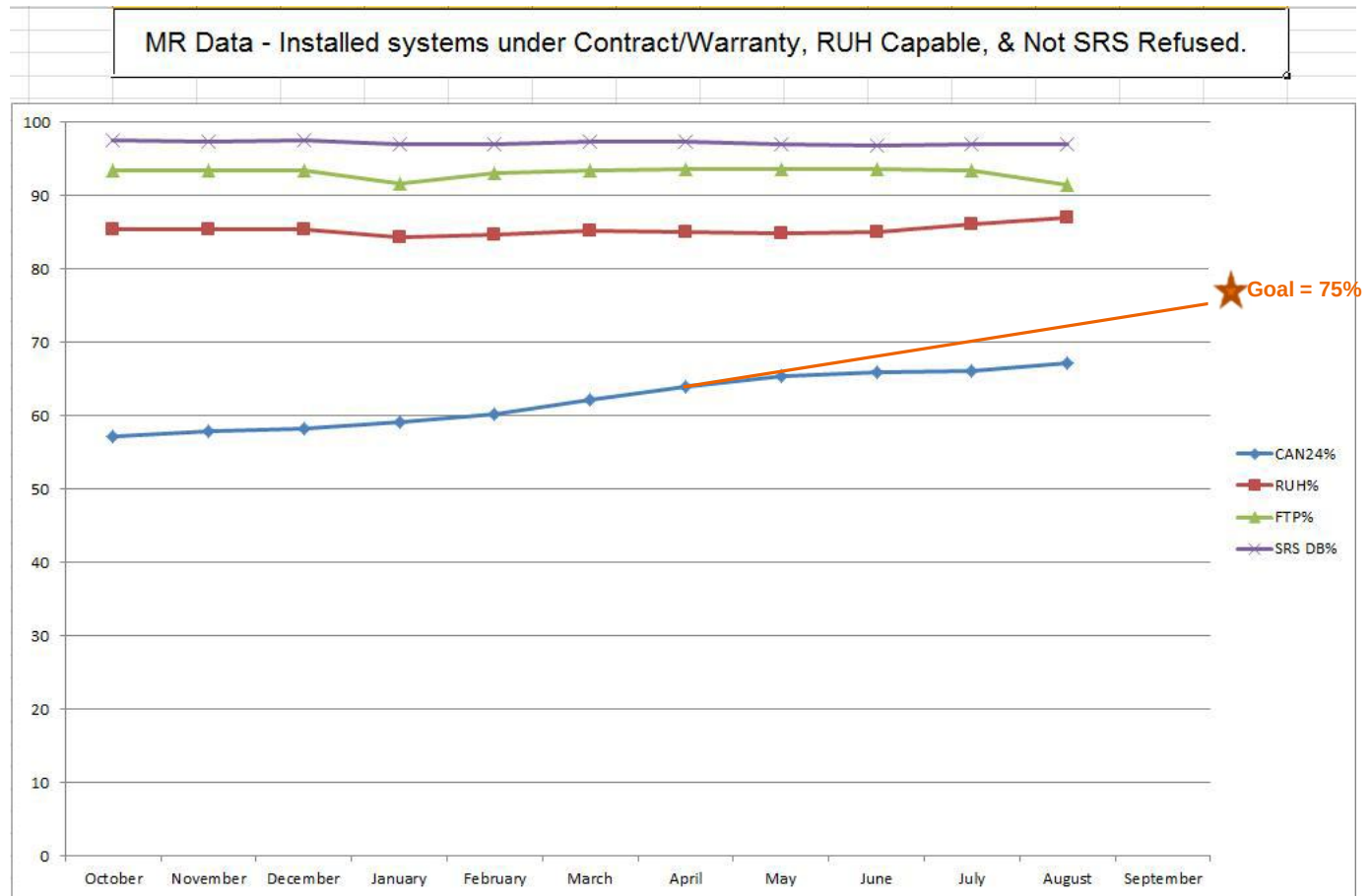

MR CAN24 Heartbeat Check to the MR Quality Page

<https://mrqprod.med.siemens.com/mrquality/mr/format/Index1.htm>



Of the 57 assigned RSM areas, 20 RSM's are currently above the 75% goal:

Here are the details:

Month	Total CAN24's	CAN24's Sending Heartbeat Check	%
October	--	--	--
November	1289	746	57.9%
December	1319	767	58.2%
January	1339	789	58.9%
February	1376	831	60.4%
March	1398	869	62.2%
April	1420	907	63.9%
May	1358	889	65.5%
June	1384	913	66.0%
July	1417	937	66.1%
August	1441	969	67.2%
September			
GOAL----->			75.0%

SRS MR_CAN24 Heartbeat Check Report

Top Performing RSM's


Rank	ZVP	RSM	Zone	SRS Rep	CAN24's	CAN24's FTP?	%
1	Campbell, John	Christian Maldonado, Jose	SOUTHEAST ZONE	Valdez Ladd	1	1	100.0%
2	Campbell, John	Green, Andrew	SOUTHEAST ZONE	Sammy Parker	24	22	91.7%
3	Parrella, Mike	Pedersen, Dale	CENTRAL ZONE	Ron Anderson/Valdez Ladd	37	33	89.2%
4	Lawrie, Stuart	Open RSM055	WESTERN ZONE	Konrad Grygoruk	9	8	88.9%
5	Campbell, John	Bruner, Perry	SOUTHEAST ZONE	John Acrivides	16	14	87.5%
6	Campbell, John	Fultz, David	SOUTHEAST ZONE	Sammy Parker	15	13	86.7%
7	Parrella, Mike	Rock, James	CENTRAL ZONE	Rick Hendon	21	18	85.7%
8	Pasquin, Victor	Caron, Glen	NORTHEAST ZONE	William Knight	33	28	84.8%
9	Pasquin, Victor	Irshad, Faisal	NORTHEAST ZONE	William Knight	17	14	82.4%
10	Pasquin, Victor	Casey, Sean	NORTHEAST ZONE	William Knight	22	18	81.8%
11	Parrella, Mike	Chadwick, John	CENTRAL ZONE	AJ Demianenko	16	13	81.3%
12	Campbell, John	Schwalje, Michael	SOUTHEAST ZONE	Sammy Parker	30	24	80.0%
12	Parrella, Mike	Kleve, Ronald	CENTRAL ZONE	Valdez Ladd	15	12	80.0%
12	Campbell, John	Saunders III, Robert	SOUTHEAST ZONE	Sammy Parker	15	12	80.0%
15	Parrella, Mike	Jeffery, Jeremy	CENTRAL ZONE	Valdez Ladd	52	41	78.8%
16	Campbell, John	Wilson, Randolph	SOUTHEAST ZONE	John Acrivides	46	36	78.3%
17	Lawrie, Stuart	Nasrati, Randy	WESTERN ZONE	Ron Anderson/Konrad Grygoruk	59	46	78.0%
18	Werkmeister, Kenneth	Sontag, Douglas	MIDWEST ZONE	Sam Moore	17	13	76.5%
19	Parrella, Mike	Weber, Roy	CENTRAL ZONE	AJ Demianenko/Sam Moore	25	19	76.0%
20	Parrella, Mike	Palmer, Michael	CENTRAL ZONE	Valdez Ladd	29	22	75.9%

MR New Installs with CAN24 tested BEFORE the turnover date

YEARLY TOTALS									
MR_CAN24 Heartbeat Check BEFORE Turnover Date									
***** COMPLETELY UPDATED *****									
Month	Total Systems	CAN24 Tested	CAN24 %		Zone	Total Systems	CAN24 Tested	CAN24 %	
OCTOBER	10	6	60.0%		W	41	21	51.2%	
NOVEMBER	25	4	16.0%		NE	43	20	46.5%	
DECEMBER	23	7	30.4%		C	64	28	43.8%	
JANUARY	17	6	35.3%		MW	72	30	41.7%	
FEBRUARY	30	18	60.0%		S	60	21	35.0%	
MARCH	39	23	59.0%		Totals	280	120	42.9%	
APRIL	18	10	55.6%						
MAY	27	14	51.9%						
JUNE	27	8	29.6%						
JULY	31	10	32.3%						
AUGUST	33	14	42.4%						
SEPTEMBER									
Totals	280	120	42.9%						
					Reason CAN24 Not Tested at time of Install	Total Missed	Pct		
					1 = New Construction	35	21.9%		
					4 = ISE did not test CAN24	32	20.0%		
					2 = Network drop not in room	28	17.5%		
					5 = Customer Firewall/Routing/Ports Issue	28	17.5%		
					3 = Network drop is available, but no IP requested	19	11.9%		
					7 = CAN24 Technical Issue	10	6.3%		
					6 = SRS Refused/DOD	8	5.0%		
					Totals	160			

SRS CONNECTION TSE STATE ASSIGNMENT MAP

Version - 21
2018

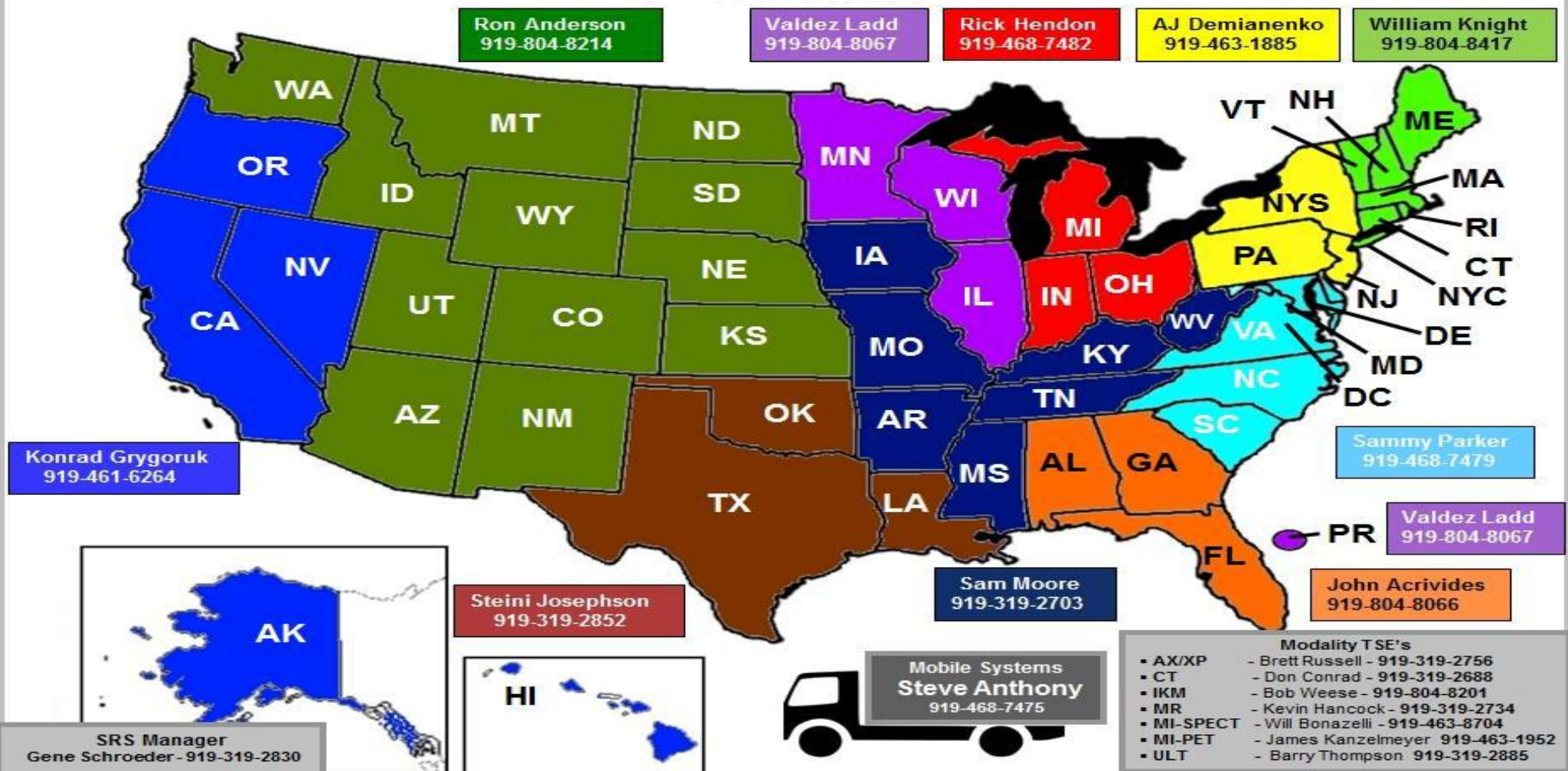

Agenda:

- Safety First!
- MR – Quality Initiatives
- Updates
- Speed Infos
- Tips of the month
- Max Update
- Special Topic
- **Feedback from field**

Questions?