

MAGNETOM

Software System

Software Installation (VA25A)

Document Version

Siemens reserves the right to change its products and services at any time.

In addition, manuals are subject to change without notice. The hardcopy documents correspond to the version at the time of system delivery and/or printout. Versions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order a current version or refer to our website <http://www.healthcare.siemens.com>.

Disclaimer

Siemens provides this documentation “as is” without the assumption of any liability under any theory of law.

The content described herein requires superior understanding of our equipment and may only be performed by qualified personnel who are specially trained for such installation and/or service.

Copyright

“© Siemens, 2004” refers to the copyright of a Siemens entity such as:

Siemens Healthcare GmbH - Germany

Siemens Aktiengesellschaft - Germany

Siemens Shenzhen Magnetic Resonance Ltd. - China

Siemens Shanghai Medical Equipment Ltd. - China

Siemens Healthcare Private Ltd. - India

Siemens Medical Solutions USA Inc. - USA

Siemens Healthcare Diagnostics Inc. - USA and/or

Siemens Healthcare Diagnostics Products GmbH - Germany

1	General	7
1.1	Introduction.	7
1.2	syngo MR system software	8
1.2.1	Prerequisites	8
1.2.2	MR 2004A Rev. B	8
1.2.3	How to use this Register	8
1.2.4	Replacing the HOST.	8
1.3	CCU - Firmware / Loadware / Testware	10
1.3.1	Cameras.	10
1.4	Customer files	11
1.5	Remote diagnostics	12
1.6	Safety	13
1.6.1	SAR	13
1.6.2	Stimulation Monitor	13
1.6.3	Access to the system	13
1.6.4	Integrity of software back-up	13
1.7	How to Navigate in this Manual	15
1.7.1	Start points	15
1.7.2	DICOM Configuration	16
2	Getting started and configuration pages	17
2.1	General	17
2.1.1	Passwords	17
2.1.2	Prerequisites for the main console	17
2.1.3	Prerequisites for the satellite console	17
2.2	Site specific configuration.	18
2.2.1	Prerequisites MRC	18
2.2.2	Prerequisites MSC	18
2.3	Configuration pages for MRC	20
2.3.1	System.	20
2.3.2	Database	24
2.3.3	Service.	24
2.3.4	DICOM	28
2.3.5	Import/Export	34
2.3.6	EPR Interface	35
2.3.7	External Devices	35
2.3.8	Applications.	35
2.3.9	Measurement	38
2.4	Configuration pages MRSC	40
2.4.1	System.	40
2.4.2	Local host	41
2.4.3	Database	43
2.4.4	DICOM	45
2.4.5	Import/Export	45
2.4.6	External Devices	47
2.4.7	Applications.	47

3	Installing MRC / MRC+MRSC	51
3.1	Overview of installing VA25A MRC/MRC+MRSC	51
3.2	General	53
3.3	Installing Software	54
3.3.1	Prerequisites	54
3.3.2	Installing MR 2004A Rev. B MRC / MRSC	55
3.3.3	Keyboard settings and Digit grouping symbol for administrator	56
3.3.4	Export MRC configuration data to MRSC	60
3.3.5	Import configuration data from CD-R to MRSC	61
3.3.6	Reboot	62
3.4	Initial configuration after software installation MRC and MRSC	63
3.4.1	Keyboard settings and Digit grouping symbol for meduser	63
3.5	MRC	68
3.5.1	Activate BIOS protection	68
3.6	IFN (Image Filter New) Installation Work Steps	69
3.6.1	Hints	69
3.6.2	Installing IFN on MRC and MRSC	70
3.6.3	Verifying correct IFN installation	71
3.6.4	CV License (only valid for systems delivered before 01/2006)	71
3.7	MRSC	72
3.7.1	General	72
3.7.2	MRSC Sharing Database with MRC	73
3.7.3	Completed configuration	73
3.7.4	Activate BIOS protection	73
3.8	IFN (Image Filter New) Installation Work Steps	74
3.8.1	Hints	74
3.8.2	Installing IFN on MRC and MRSC	75
3.8.3	Verifying correct IFN installation	76
3.8.4	CV License (only valid for systems delivered before 01/2006)	76
4	Installing Software at the Imager	77
4.1	Overview installing Imager SW VA25	77
4.1.1	General	77
4.2	Installing Software at the Imager "R630"	78
4.2.1	Prerequisites	78
4.2.2	Installation procedure	78
5	Site specific data	80
5.1	Save site specific data	80
5.1.1	Numaris4	80
5.1.2	CustomerProtocols	80
5.1.3	SW-Settings	80
5.1.4	Prerequisites	80
5.1.5	Settings that are not part of the backup/restore	80
5.1.6	Backup (Save site specific data) on CD	81
5.1.7	Backup (Save site specific data) on USB device	82

5.1.8	Additional packages	83
5.2	Restoring site specific data MR 2004A Rev. B on MR 2004A Rev. B	84
5.2.1	Prerequisites	84
5.2.2	Restore site-specific from CD	84
5.2.3	Restore site-specific from external USB device	86
5.2.4	Settings that must be manually activated after restore	86
6	Configuration	89
6.1	General	89
6.2	Prerequisites	90
6.3	Access Service Software	91
6.3.1	Changes in the configuration	91
6.4	Parallel Printer Installation	92
6.4.1	Installing local parallel printer	92
6.4.2	Enabling MRC to access printer share of MRSC	92
6.5	Network Printer	93
6.5.1	Installing network printer	93
6.6	General Printer Properties	94
7	Remote connectivity	99
7.1	General	99
7.2	Prerequisites	100
7.2.1	System Management	100
7.2.2	Pro Active Service	100
7.2.3	Quality Assurance DB (QDB)	100
7.2.4	Event DB (EDB)	100
7.2.5	System Utilization	100
7.2.6	Check / add route to SIEMENS Remote Services router	100
7.3	Configure MAIL and FTP target (router systems)	102
7.4	Testing connectivity from system to Rdiag server	104
7.4.1	Send a test mail	104
7.4.2	Test FTP login to the Remote Server	105
7.5	AutoReport configuration	107
7.5.1	Configuring QDB report	107
7.5.2	Configuring EDB report	108
7.5.3	Configure System Utilization reporting	109
7.6	System Utilization	111
7.7	System Management Installation	112
7.7.1	Prerequisites	112
7.7.2	Routing	113
7.7.3	Check connectivity to the System Management server	114
7.7.4	Required ports for System Management (CA Unicenter) Com- munication	114
7.7.5	Master	114
7.8	Troubleshooting and useful information	117
7.8.1	System Management Installation failed	117

8	Appendix	120
8.1	Re installation of software fails	120
8.1.1	Formatting the system disk	120
8.2	System language and weight	122
8.3	Keyboard	123
8.4	Language Settings Japanese	124
8.4.1	Switch language from international (english-US) to japanese . . .	124
8.4.2	Switch language from japanese to international (english-US) . . .	127
8.5	Language Settings Chinese	130
8.5.1	Switch language from international (english-US) to chinese . . .	130
8.5.2	Switch language from chinese to international (english-US) . . .	133
8.6	Institution name	135
8.7	Regional settings	136
8.7.1	Date and Time syntax	136
8.8	System Manager	137
8.8.1	Restart Application	137
8.9	Copy SIEMENS Programs to the customer tree	138
8.10	IP Security	139
8.10.1	Adding host to IP security list	139
8.10.2	Deleting all IP security settings	139
8.10.3	Adding all configured hosts to IP security list	139
8.11	Context Vision	141
8.11.1	Requested data for cvlicense generation	141
8.11.2	Activate license for Context Vision	141
8.12	MR auto login	143
8.12.1	Access NT login	143
8.13	syngo startup / shutdown counter	144
8.13.1	How to terminate the syngo counter	144
8.14	Terminate deadlock of crashed Service UI	145
8.14.1	How to terminate the Service IU (Internet Explorer).	145
8.15	Expert - i	146
8.16	Measurement system configuration matrix	147
8.17	Dicom	153
8.17.1	Dicom and UNIX nodes	153
8.17.2	Dicom configuration data for SIEMENS products	153
9	Changes to previous version	156
10	List of Hazard IDs	157

1.1 Introduction

The installation and service of equipment described herein is to be performed by qualified personnel which is employed or authorized by Siemens Medical Solutions or one of its affiliates.

Assemblers and other personnel not employed by nor directly affiliated with nor authorized by Siemens Medical Solutions or one of its affiliates technical services are directed to contact the local office of Siemens Medical Solutions or its affiliates before attempting installation or service procedures.

1.2 syngo MR system software

1.2.1 Prerequisites

To configure the syngo-based MR system following training is required in combination with this Software manual:

- XP-administrator training
- Network training (TCP/IP)
- DICOM training
- syngo training / syngo MR training

1.2.2 MR 2004A Rev. B

The following software media must be delivered with the system: **CD for VA25A**

The MR 2004A Rev. B software includes the service software which is protected by a license mechanism.



The "Release Notes syngo MR MR 2004A Rev. B" provides additional information regarding the MR 2004A Rev. B software.

1.2.2.1 MRC / MRSC R630

Operating system Windows XP Professional (syngo MR 2004A Rev. B)

1.2.2.2 Imager R630

Operating system Windows XP Professional (syngo MR 2004A Rev. B)

1.2.3 How to use this Register

Save time by reading this chapter first

This chapter contains reference, safety, and information regarding navigation in this manual.

1.2.4 Replacing the HOST

- The dongle from the defective host has to be installed at the new host due to the license mechanism.

- Applications licenses are generated in relation to the flexid of the 7, 9- (Dongle) with Census plus.

Exception: if a 8- Dongle is plugged, all application licenses are generated in relation to the MAC address of the first NIC (one exception is the Context Vision filter, it needs two licenses).

The low level license of CV filter is always generated in relation to the DONGLE.

The Service license is generated in relation to the DONGLE always (7, 9- and 8- dongles).

- IP settings from the own host might not be restored from the save set, if the old NIC PCB differs to the new NIC: Enter the local IP address manually in the service platform.
- If the Graphic card differs in between the original host and the new host: Perform the monitor calibration.



All licenses for customer options are handled by TD HSC 1 (phone +49 9191-188080).

1.3 CCU - Firmware / Loadware / Testware

1.3.1 Cameras



The use of non-released combinations of hardcopy cameras and BU products is prohibited and does not meet the requirements defined. For detailed information regarding the CAMERA configuration, refer to technical documentation "Camera Info Numaris 4 VA25A"

1.4 Customer files

User settings for Filming, Viewing, etc., are part of the save set package "**User-settings**" in SW-Settings_02.

The Customer Protocols are stored in the save set package "**CustomerProtocols**" (Customer protocols are part of the Numaris4 save set package as well).

Ask the customer to contact service for support after making major changes to user-defined protocols and user settings.

1.5 Remote diagnostics

Only IE 6.0 has been tested and released for remote service as service browser at syngo MR systems running with NUMARIS MR 2004A Rev. B.

The "SIEMENS Remote Services TM" will only support router connections for NEWLY delivered systems. See (→ Siemens Remote Service, <http://cs.med.siemens.de/cms/productinfo/default.asp?startfolder=e3371dfd0ace1e5181c7d36777ef3b4b>) for the most recent information regarding the router configuration..

(H WS CS - Customer Services -> For Service -> Product Information -> Siemens Remote Service)

The project manager is responsible to plan for a router connection to the Siemens Remote Diagnostic server. Regarding questions to your pre configured router from Siemens please contact mailto:srs_final@med.siemens.de or your remote diagnostic responsible.

1.6 Safety

1.6.1 SAR

The absorption of RF-energy by the patient is an integral part of every MR examination. The amount of energy absorbed has to be within the limits specified by the respective federal agency.

The Specific Absorption Rate (SAR) is different for each country. For this reason, you have to select the country code in the site-specific configuration corresponding to the local guidelines and regulations.

In addition, the temperature and relative humidity in the examination room affect the SAR values (only applies to countries subject to FDA and IEC guidelines, that is, all countries except Japan).

1.6.2 Stimulation Monitor

The magnetic field changes very quickly on Sonata and Quantum systems during gradient switching. If the electric field generated exceeds a specific threshold value, an electric current may be induced in the patient.

The Stimulation Monitor threshold is different for each country. For this reason, you have to select the country code in the site-specific configuration corresponding to the local guidelines and regulations.

The current induced may cause noticeable peripheral nerve stimulations and a feeling of discomfort for the patient.

1.6.3 Access to the system

The system software can only be accessed via password-protected accounts. This safety feature protects confidential (patient) data and reduces the possibility of (accidental) software manipulation by a third party.

The customer is provided with the passwords for the accounts required to operate and maintain the system.

1.6.4 Integrity of software back-up

The software has undergone extensive testing. However, this does not ensure that the software will run indefinitely without errors.

A backup (save set of site-specific data) should be created to avoid loss of system-specific data if the software has to be restored, e.g. after a system disk crash.

The save set will reduce system downtime to a minimum.



Before you release the system to the customer, ensure that the image quality is within specifications and that all safety-relevant checks have been passed, e.g., the image orientation is correct.

Refer to Register "Tune-up" in the Technical Documents.

1.7 How to Navigate in this Manual

This manual includes information about various tasks. Please see the corresponding chapters listed in the table:

1.7.1 Start points

Getting started and configuration pages Your notes, workflow for initial configuration of new systems (software is pre-installed at the factory)	Refer to (→ Getting started and configuration pages / Page 17) 
Installing MRC / MRC+MRSC MRC or MRSC suffers from a defect that requires re-installation of the "syngo MR 2004A Rev. B"	Refer to (→ Installing MRC / MRC+MRSC / Page 51) 
	Refer to (→ Installing MRC / MRC+MRSC / Page 51) 
Installing Software at the Imager The Imager suffers from a defect that requires re-installation of the "syngo VA25A"	Refer to (→ Installing Software at the Imager / Page 77) 
Site specific data Save or Restore site-specific data, e.g. to save the configuration after making changes or restore site-specific data after software installation	Refer to (→ Site specific data / Page 80) 
Configuration Configuring / reconfigure a running system.	Refer to (→ Configuration / Page 89) 

Remote connectivity Configuring Quality Assurance Reporting and System Utilization reporting.	Refer to (→ Remote connectivity / Page 99) 
Appendix General information (e.g., BIOS settings and CV license, Institution name, Date, Time, Measurement system data.	

1.7.2 DICOM Configuration

The "CCC Team" provides a powerful configuration management tool see <http://cs.med.siemens.de/> "For Service", "DICOM", in the DICOM Homepage (CCC) navigate to "Interoperability DB", "Download", "Acquisition Tool".

2.1 General

This chapter describes the steps involved in configuring / reinstalling syngo-MR systems running VA25A.

MRC/MRSC are used as abbreviations for the Main/Satellite consoles, mrc/mrsc are used as host names.

This part of the table is always valid.
This part of the table is only valid for MRSC consoles

2.1.1 Passwords



All passwords required to install and configure the system are available on the Med TD Intranet server at <http://www-td.med.siemens.de/>, "CS Knowledge Base", Speed Infos, search in "Speed Info" (except for the service key = variable password, which is handled by the districts).

- For more information regarding the service key (variable password), refer to the "Service Software" in TSG.

2.1.2 Prerequisites for the main console

- The VA25A software is installed at the MRC.

2.1.3 Prerequisites for the satellite console

- The VA25A software is installed at the MRC and MRSC.



If all prerequisites are fulfilled, choose one of the suggested procedures and follow the instructions and the flowchart.

2.2 Site specific configuration

Recommendation: For paper versions, take out this page and use it for navigation.

2.2.1 Prerequisites MRC

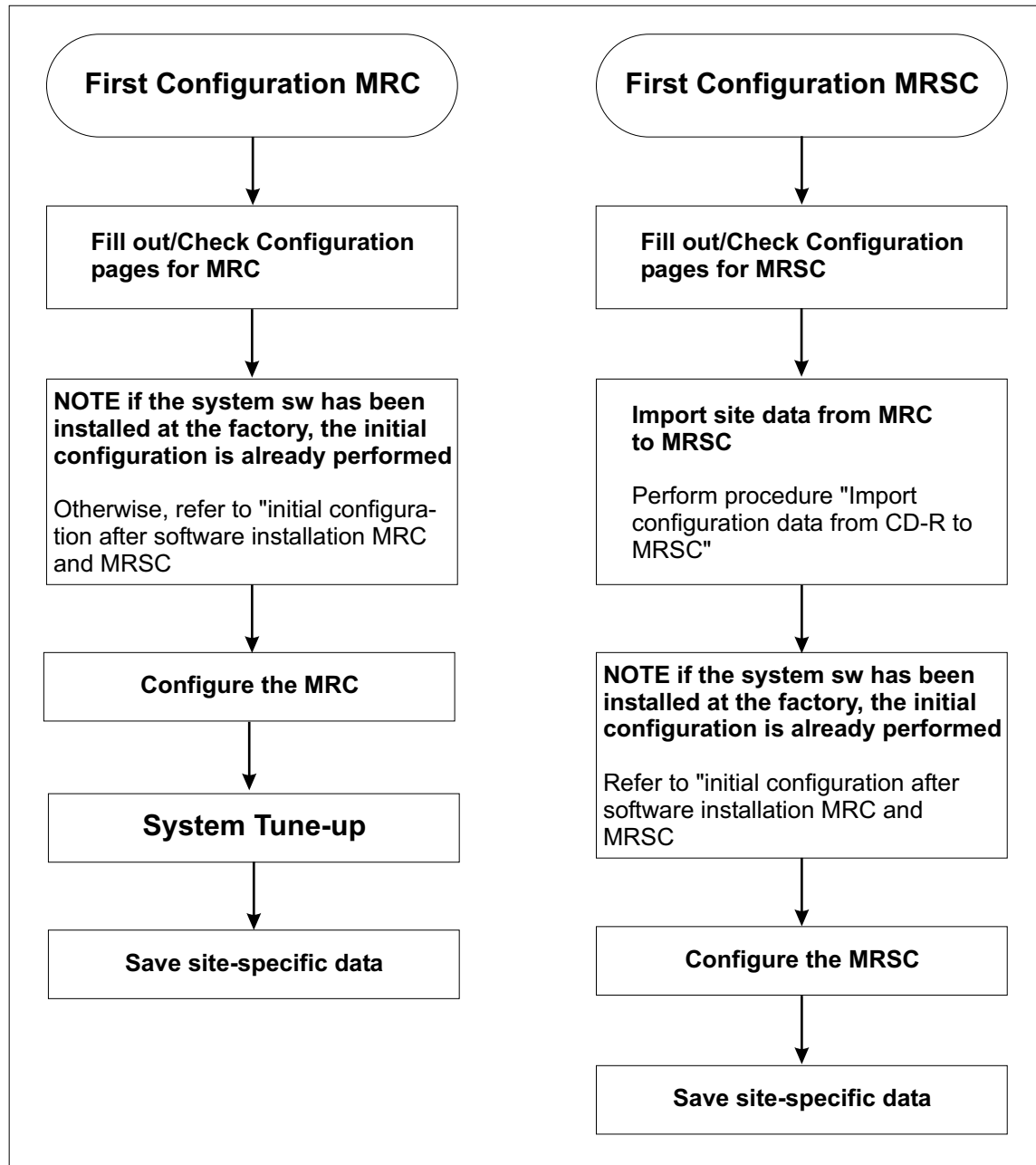
- System software MR 2004A Rev. B is installed at the MRC

2.2.2 Prerequisites MSC

- Configuration is completed for the MRC

- System software MR 2004A Rev. B is installed at the MRSC

Fig. 1: First Configuration MRC/MRSC



2.3 Configuration pages for MRC

Before you start the configuration of the system, please fill out the following configuration pages. The order of the tables in the configuration pages corresponds to the order during "Full Configuration".

Actual SW version	
Software version	Signature / Date

2.3.1 System



If the system is connected to a hospital network, contact the network administrator regarding IP addresses and network settings.

list of system options	
Option / System type	Your notes
Attached to network	
Modem	
Main Console with Satellites or Satellite Console	
DICOM Print Devices	
DICOM Offline Devices	
DICOM Networking	
DICOM HIS/RIS	
Image Import/Export	
Paper Printer	
EPRI	

2.3.1.1 Local Host

MR System Specification		
Input field	Your notes	Note
Console Type / System Type		
System type e.g.: Harmony, Symphony, Sonata, Concerto, HarmonyExpert, SymphonyVision, SonataVision, Allegra, or Trio		

Site Info		
Input field	Your notes	Note
System		
Serial No.		
Maintenance Unit	mr	
Department		
Service		
Site Identification No		
Service Center Phone		
Customer and Address		
Name		
Id		
Hospital		
Street		
Street No		
Zip Code		
Phone No		
City		
District		
Country		
Time zone		
		Pull-down menu

Site Info		
Input field	Your notes	Note
AutoDaylightTime		
(Day light saving)		Check box



The host name must be unique in the network!

TCP/IP LAN		
Input field	Your notes	Note
Identification		
Computer name	Recommended syntax: e.g. MRC22000 or MRSC22000	(MRC and system serial number)
		Host name for MRC Use only letters from a to z (numbers are allowed as well, but they may not be used for the first character of the host name)
Adapter (Network adapter)		
select name		Pull-down menu lets you select from a predefined list
IP Address		
IP address (local host)		Delete the defaults in the configuration menu for the first Network PCB (only one entry is supported in this version)
Subnet mask		
gateways		
NOTE: Use default addresses for the second Network PCB 192.168.2.1 (NIC for measurement system network MPCU and Imager)		
WINS		
Primary WINS server		
Secondary WINS server		
Enable DNS for Windows Resolution		Select box

TCP/IP LAN		
Input field	Your notes	Note
Enable LMHOSTS Lookup		Select box
DNS		
Domain		
1 DNS service search order		Several DNS nodes can be defined
2 DNS service search order		
3 DNS service search order		

Users		
Input field	Your notes	Note
Account for UI	Leave default Default is meduser	Changing the password / account will result in problems during manual adjustments
Default Password / or give	Leave default password	
Account for PPP (Account name is: PPP)		
Default Password / or give		Account name: PPP This account is required for remote service sessions via modem. You must agree on a unique PPP password for your system together with the responsible USC/HSC Rdiag server administrator
Enable Autologin		Autostart of the MR application

Monitor Type settings		
Input field	Your notes	Note
Monitor Type		
color / greyscale		For color LCD monitors select: Appearance Siemens (True Color)

2.3.2 Database

Satellite Console (host name and IP-address of the MRSC)		
Input field	Your notes	Note
Host Properties (MRSC)		
Select Host		Enter the IP address of the MRSC and host name Only LAN is supported!

2.3.3 Service

Autoreport (mail target settings for service autoreporting tool)		
Input field	Your notes	Note
setting of LAN connection of remote mail server		
Select Host		Several target hosts can be set Note the logical name
Host Properties (remote server host)		
1 Host Name		Note the host names and IP addresses
1 TCP/IP address		
2 Host Name		
2 TCP/IP address		
3 Host Name		
3 TCP/IP address		
4 Host Name		
4 TCP/IP address		

Mail		
Input field	Your notes	Note
Host Properties		
Select Host		Enter the IP address of the Rdiag server

Mail		
Input field	Your notes	Note
1 Host Name		Note the host names and IP addresses
1 TCP/IP address		
2 Host Name		
2 TCP/IP address		
Mail Server Host Selection		
email address		Enter the E-mail address of your local USC
email address in Remote Service Access Control window		Check-box (If this check-box is set the customer function Mail in the "Remote Access Control window" will be sent to the above e-mail address)

FTP (ftp target settings for service file&image transfer tool)		
Input field	Your notes	Note
Select Host		Several target hosts can be set Note the logical name
Host Properties (Remote server host)		
1 Host Name		Note the host names and IP addresses
1 TCP/IP address		
2 Host Name		
2 TCP/IP address		
Logical Name		
FTP Host Selection		Pull-down menu
FTP User		
FTP Password		(Confirm password at the second row)

AutoTransfers		
Input field	Your notes	Note
activate/deactivate automatic Transfers (initiated by the autoreport tool)		
activate automatic Transfers		Must be enabled that Auto Reports can be sent automatically by the system. system utilization will only work if check box is enabled.

Eventlog		
Input field	Your notes	Note
Configuration of display parameters for the Eventlog		
Log only Service messages		Check box
Display Service messages		Check box
Display sequence number		Check box

Backup/Restore		
Input field	Your notes	Note
package		
Additional package 1		The following packages are pre selectable by default: Numaris4 SW-Settings CustomerProtocols Diagnostic ExamDB Additional packages might be included for re
Additional package 2		
Additional package 3		
Additional package 4		
Additional package 5		
Additional package 6		
groups		

Backup/Restore		
Input field	Your notes	Note
files & directories		
file & dir		
backup pre		
backup post		
restore pre		
restore post		

Licensing		
Input field	Your notes	Note
Input/Change of Licenses		
License Key		Keep the list with licences separated from this documentation
Hardware ID Type	Dongle	
Hardware ID		Unique dongle number (ID)

Fonts		
Input field	Your notes	Note
Input/Change of Licenses		

Fonts		
Input field	Your notes	Note

2.3.4 DICOM

The "CCC Team" provides a powerful configuration management tool (see "DICOM", CCC Connecting Systems, "Databases", "Interoperability db", "Download", "Acquisition Tool")

(→ Acquisition Tool, https://intranet.healthcare.siemens.com/cms/IKM/en/Med-wideTopics/Connectivity_Competence_Center/Databases/Interoperability_DB/Pages/Acquisition_Tool.aspx)



Ask your hospital network administrator for the correct DICOM nodes.

General		
Input field	Your notes	Note
Local application entity titles		
HIS / RIS		(SCU)
Study Transfer		(SCP + SCU)
Print		(SCU)
Simultaneous DICOM associations		
maximum PDU size		(Packet Data Unit size)
time-out values	use default settings	
implicit raw data transfer		Check-box

character set (configuration of the DICOM default character set used when registering new patients)

Input field	Your notes	Note
DICOM default character set		

character set (configuration of the DICOM default character set used when registering new patients)

Input field	Your notes	Note
	Set radio button to "allow multiple character sets"	ISO 2022 IR100 Japan Only: To support Japanese character set, select ISO 2022 IR13
Select character set		Pull-down menu

Offline devices (Setting of DICOM archive devices on local host)

Input field	Your notes	Note
Select Logical Name		
Logical name		Several targets can be set e.g CDRW MOD etc. For MOD drives: set the SPI file system flag
Logical name		
Logical name		
Logical name		

Network Nodes setting of LAN connection of remote DICOM transfer target)

Input field	Your notes	Note
Host Properties		
Select Host	1	Several target hosts can be set Note the logical name
	2	
	3	
	4	
	5	
	6	
Host Properties (remote server host)		

Network Nodes setting of LAN connection of remote DICOM transfer target)		
Input field	Your notes	Note
Host Name	TCP/IP address	Note the host names and IP addresses
1		
2		
3		
4		
5		
6		

Network Nodes (definition of application entity titles and ports of DICOM storage SCPs)		
		Note
General Node Properties		
Logical Name	(Default node / syngo based)	(syngo based yes/no, first default node or second default node etc.)
1.	1.	
2.	2.	
3.	3.	
4.	4.	
5.	5.	
6.	6.	
Application Entity Properties		
AE Title	Port Number	
1.	1.	
2.	2.	
3.	3.	
4.	4.	
5.	5.	
6.	6.	
Supported DICOM services		
Storage, Storage, Commitment, Query, and/or Retrieve		

Network Nodes (definition of application entity titles and ports of DICOM storage SCPs)		
		Note
1.		If you are unsure which services are needed, please select all
2.		
3.		
4.		
5.		
6.		

Print Devices (definition of application entity titles and ports of DICOM storage SCPs)		
		Note
General Node Properties		
Logical Name	(Default node / syngo based)	(syngo based yes/no, first default node or second default node etc.)
1.	1.	
2.	2.	
3.	3.	
4.	4.	
5.	5.	
6.	6.	
Application Entity Properties		
AE Title	Port Number	
1.	1.	
2.	2.	
3.	3.	
4.	4.	
5.	5.	
6.	6.	
Supported DICOM services		
Print		

hc overview (overview of hardcopy devices)		
Input field	Your notes	Note
Select logical name 1		
Name		Remember to mark the Default Camera
Type		
Comment		
Location		
Hold printed filmjobs		(Pulldown menu)
Network access		
DICOM access		
Select logical name 2		
Name		
Type		
Comment		
Location		
Hold printed filmjobs		(Pulldown menu)
Network access		
DICOM access		
Select logical name 3		
Name		
Type		
Comment		
Location		
Hold printed filmjobs		(Pulldown menu)
Network access		
DICOM access		



The import of LUT files is only required in certain cases.

Image type / Max Density		
Input field	Your notes	Note
LUT	MR1	Use defaults
Density	3.00	(NEW with MR 2004A Rev. B LUT with 16 bit are available Presentation LUT is supported.

HIS/RIS nodes (setting of LAN connection to remote HIS/RIS station)		
Input field	Your notes	Note
Select Host		Several target hosts can be set Note the logical name
Host Properties (remote server host)		
1 Host Name		Note the host names and IP addresses
1 TCP/IP address		
2 Host Name		
2 TCP/IP address		
3 Host Name		
3 TCP/IP address		
4 Host Name		
4 TCP/IP address		
Select Logical Name		
Logical name		Several targets can be set e.g. CDRW MOD etc. For MOD drives: set the SPI file system flag
Logical name		
Logical name		
Logical name		
definition of application entity title and port of DICOM report or worklist SCP		
General Node Properties		
Logical Name		
Host		

HIS/RIS nodes (setting of LAN connection to remote HIS/RIS station)		
Input field	Your notes	Note
Application Entity Properties		
AE Title		
Port Number		
Supported DICOM services		
Basic Worklist		Check box
Performed Procedure Step		Check box

2.3.5 Import/Export

directories (setting of LAN connection to remote computer)		
Input field	Your notes	Note
Host Properties		
Select Host	1	Several target hosts can be set Note the logical name
	2	
	3	
Host Properties (remote server host)		
Host Name	TCP/IP address	Note the host names and IP addresses
1		
2		
3		
Select logical HD device		
Logical Name	Path	Drive
1		
2		
3		
Remote system		

directories (setting of LAN connection to remote computer)		
Input field	Your notes	Note
Host	Share	
1		
2		
3		

2.3.6 EPR Interface

Server		
Input field	Your notes	Note
Location of EPRServer		http://123.456.789.10/ webEPR
Type of EPRServer	SHS	

2.3.7 External Devices

Paper Printer		
Input field / Your notes		
Patient		
Letter small:	A3:	A4:
A3:	B4:	B5:
A5:	Legal:	

2.3.8 Applications

Pat. Registration		
Input field		Your notes
Displayed Patient Registration Dialog Elements		
Patient	Study	
Name Fields	Study List	
Age	Accession No	
Height	Requested Procedure	

Pat. Registration		
Input field		Your notes
Weight	Request Id	
Address	Study Comment	
Other Name(s)		
Other ID(s)		
Ethnic Group		
Military Rank		
Details	Visit	
Additional Info	Referring Physician	
Institution	Requesting Physician	
	Admitting Diagnosis	
	Ward	
	Admission Id	
Patient Registration Dialog Invocation		
implicitly in scheduler		

Viewer	
Input field	Your notes
Viewer configuration settings	
Close all studies when new one is loaded	
Close all series when new one is loaded	

CorRea	
Input field	Your notes
Displayed CorRea Dialog Elements	

CorRea		
Input field		Your notes
Patient	Instance	
Name Fields	Number	
Other Name(s)	Date	
Other ID(s)	Time	
Ethnic Group	Contrast Medium	
Military Rank	Comment	
Height	Calibration Image	
Weight	Image Laterality	
Address	Series	
Additional Info	Request ID	
Hospital	Modality	
Referring Physician	Patient Position	
Requesting Physician	Requesting Physician	
Admitting Diagnose	Institution Address	
Ward	Performing Physician	
Admission ID	Operator	
Study	Body Part Examined	
Accession Number	Number	
ID	Description	
Description	Product Name	
Comments	Date	
Date	Time	
Time	Laterality	

MPPS		
Input field		Your notes
Displayed MPPS Dialog Elements		
Study	Study ID	
Study Date	Study Time	

MPPS		
Input field		Your notes
Accession No.	Description	
Modality	Station	
Status	Type	
End Date	End Time	
Dose Dialog	Comments	

Worklist Results		
Input field		Your notes
Displayed HISRIS worklist result attributes		
Date of birth	Sex	
Accession Number	Admission ID	
Referring Physician	Performing Physician	
Requesting Physician	Ward	

2.3.9 Measurement

System type		
Input field	Your notes	Note
Patient		
MPCU type:		MPCU/MPCU-2/MPCU-3:
# of RX4 boards		1 or 2
# of RF receiver		1, 2 or 4
gradient amplifier		Type of gradient AMP
gradient coil		
Stimulation Monitor		
SAR guideline		
Advanced SHIM option		
PMU option		
Unit type		(stationary / mobile)

System type		
Input field	Your notes	Note
MRIR settings		
# of Harddisks		1 or 4

Only valid for MAGNETOM Concerto (007)		
Input field	Your notes	Note
Easy site	Yes / No	Sites without external interferences, e.g. railway etc. Refer to site measurement protocol
EFI	none / installed	
mains frequency (EFI)	50 Hz / 60 Hz	Mains frequency on site
railway mains frequency (EFI)	16,7 Hz or 25 Hz	Only valid if a railway is close to the site

2.4 Configuration pages MRSC

Before you start the configuration of the system, please fill out the following configuration pages. The tables in the configuration pages are listed in the same order as they are requested during "Full Configuration" in install.



For the MRSC, only the configuration of the following items is important:

- Computer name
- IP address (local host)

Actual SW version	
Software version	Signature / Date

2.4.1 System

list of system options	
Option / System type	Your notes
Attached to Network	MUST be checked
Main Console with Satellites or Satellite Console	MUST be checked
DICOM Print Devices	
DICOM Offline Devices	
DICOM Networking	
DICOM HIS/RIS	
Paper Printer	

list of system options	
Option / System type	Your notes
Paper Printer	
EPRI	

2.4.2 Local host

MR System Specification		
Input field	Your notes	Note
Console Type		
Satellite Console	Satellite Console	



If the system is connected to a hospital network, contact the network administrator regarding IP addresses and network settings.



The host name must be unique in the network!

TCP/IP LAN		
Input field	Your notes	Note
Identification		
	Recommended syntax: e.g. MRSC22000	(MRSC and system serial number)
Computer name		Host name for MRSC Use only letters from a to z (numbers are allowed as well, but they may not be used for the first character of the host name)
Adapter (network adapter)		
select name		Pull-down menu lets you select from a predefined list
IP Address		

TCP/IP LAN		
Input field	Your notes	Note
IP address (local host)		Delete defaults in the configuration menu for the Network PCB (only one entry is supported in this version)
Subnet mask		
gateways		
WINS		
Primary WINS server		
Secondary WINS server		
Enable DNS for Windows Resolution		Select box
Enable LMHOSTS Lookup		Select box
DNS		
Domain		
1 DNS service search order		Several DNS nodes are definable
2 DNS service search order		
3 DNS service search order		

Users (Configuration of user accounts which will be needed for operation)		
Input field	Your notes	Note
Account for UI		Default is meduser
Default Password / or give	Leave the default Default is meduser	Changing the password / account will result in problems during manual adjustments
Enable Autologin	Leave default password	

Monitor Type settings		
Input field	Your notes	Note
Monitor Type		
color / greyscale		For color LCD monitors select: Appearance Siemens (True Color)

2.4.3 Database

Satellite Console		
Input field	Your notes	Note
Host Properties		Enter the hostname and IP address of the MRC (Main console) Only LAN is supported
Host Name		
TCP/IP address		
Console Type		
	Satellite Console	
Database Properties		
Name of Local Main DB		
Name of Local Scheduler DB		
Database Offset of local DBs [100..65500]	101	For multiple satellite consoles please set the offset values with proper spacing e.g. 100 for the first console and 200 for the second one. The main data base keeps offset 0

Eventlog		
Input field	Your notes	Note
Configuration of display parameters for the Eventlog		
Log only Service messages		Check box
Display Service messages		Check box
Display sequence number		Check box

Licensing		
Input field	Your notes	Note
Input/Change of Licenses		
License Key		Keep the list of licenses separate from this documentation

Licensing		
Input field	Your notes	Note
Hardware ID Type	Dongle	
Hardware ID		Unique dongle number (ID)

Fonts		
Input field	Your notes	Note
Input/Change of Licenses		



Following part of the configuration for the MRSC can be comfortably imported with the procedure (→ Export MRC configuration data to MRSC / Page 60).

Following configuration items

- LocalHost\Users
- TCPIP\LAN-Settings
- Database\SateliteConsole
- Database\DBOffset
- Applications\WorklistResults
- Applications\MPPS
- Applications\CorRea
- Applications\Viewer
- Applications\PatientRegistration
- EPRI\Server
- DICOM\ImportExportDirectories
- DICOM\HISRISNodes
- DICOM\LUTFiles
- DICOM\PrintDevices
- DICOM\NetworkNodes
- DICOM\OfflineDevices
- DICOM\General
- ExternalDevices\PaperPrinter
- Service\Licensing
- Service\Eventlog

- Service\AutoTransfers
- Service\FTP
- Service\Mail
- LocalHost\MonitorType
- LocalHost\SiteInfo

2.4.4 DICOM

Refer to your notes for MRC

2.4.5 Import/Export



Use same settings as for MRC (Only in some cases the customer needs different setting in-between MRC and MRSC)

TCP/IP LAN		
Input field	Your notes	Note
Identification		
	Recommended syntax: e.g. MRSC22000	(MRSC and system serial number)
Computer name		Host name for MRSC Use only letters from a to z (numbers are allowed as well, but they may not be used for the first character of the host name)
Adapter (network adapter)		
select name		Pull-down menu lets you select from a predefined list
IP Address		
IP address (local host)		Delete defaults in the configuration menu for the Network PCB (only one entry is supported in this version)
Subnet mask		
gateways		
WINS		
Primary WINS server		

TCP/IP LAN		
Input field	Your notes	Note
Secondary WINS server		
Enable DNS for Windows Resolution		Select box
Enable LMHOSTS Lookup		Select box
DNS		
Domain		
1 DNS service search order		Several DNS nodes are definable
2 DNS service search order		
3 DNS service search order		

directories (setting of LAN connection to remote computer)		
Input field	Your notes	Note
Host Properties		
Select Host	1	Several target hosts can be set Note the logical name
	2	
	3	
Host Properties (remote server host)		
Host Name	TCP/IP address	Note the host names and IP addresses
1		
2		
3		
Select logical HD device		
Logical Name	Path	Drive
1		
2		
3		
Remote system		

directories (setting of LAN connection to remote computer)		
Input field	Your notes	Note
Host	Share	
1		
2		
3		

2.4.6 External Devices

Paper Printer		
Input field / Your notes		
Patient		
Letter small:	A3:	A4:
A3:	B4:	B5:
A5:	Legal:	

2.4.7 Applications



Use same settings as for MRC (Only in some cases the customer needs different setting in-between MRC and MRSC)

Pat. Registration		
Input field		Your notes
Displayed Patient Registration Dialog Elements		
Patient	Study	
Name Fields	Study List	
Age	Accession No	
Height	Requested Procedure	
Weight	Request Id	
Address	Study Comment	

Pat. Registration		
Input field		Your notes
Other Name(s)		
Other ID(s)		
Ethnic Group		
Military Rank		
Details	Visit	
Additional Info	Referring Physician	
Institution	Requesting Physician	
	Admitting Diagnosis	
	Ward	
	Admission Id	
Patient Registration Dialog Invocation		
implicitly in scheduler		

Viewer	
Input field	Your notes
Viewer configuration settings	
Close all studies when new one is loaded	
Close all series when new one is loaded	

CorRea	
Input field	Your notes
Displayed CorRea Dialog Elements	

CorRea		
Input field		Your notes
Patient	Instance	
Name Fields	Number	
Other Name(s)	Date	
Other ID(s)	Time	
Ethnic Group	Contrast Medium	
Military Rank	Comment	
Height	Calibration Image	
Weight	Image Laterality	
Address	Series	
Additional Info	Request ID	
Hospital	Modality	
Referring Physician	Patient Position	
Requesting Physician	Requesting Physician	
Admitting Diagnose	Institution Address	
Ward	Performing Physician	
Admission ID	Operator	
Study	Body Part Examined	
Accession Number	Number	
ID	Description	
Description	Product Name	
Comments	Date	
Date	Time	
Time	Laterality	

MPPS		
Input field		Your notes
Displayed MPPS Dialog Elements		
Study	Study ID	
Study Date	Study Time	

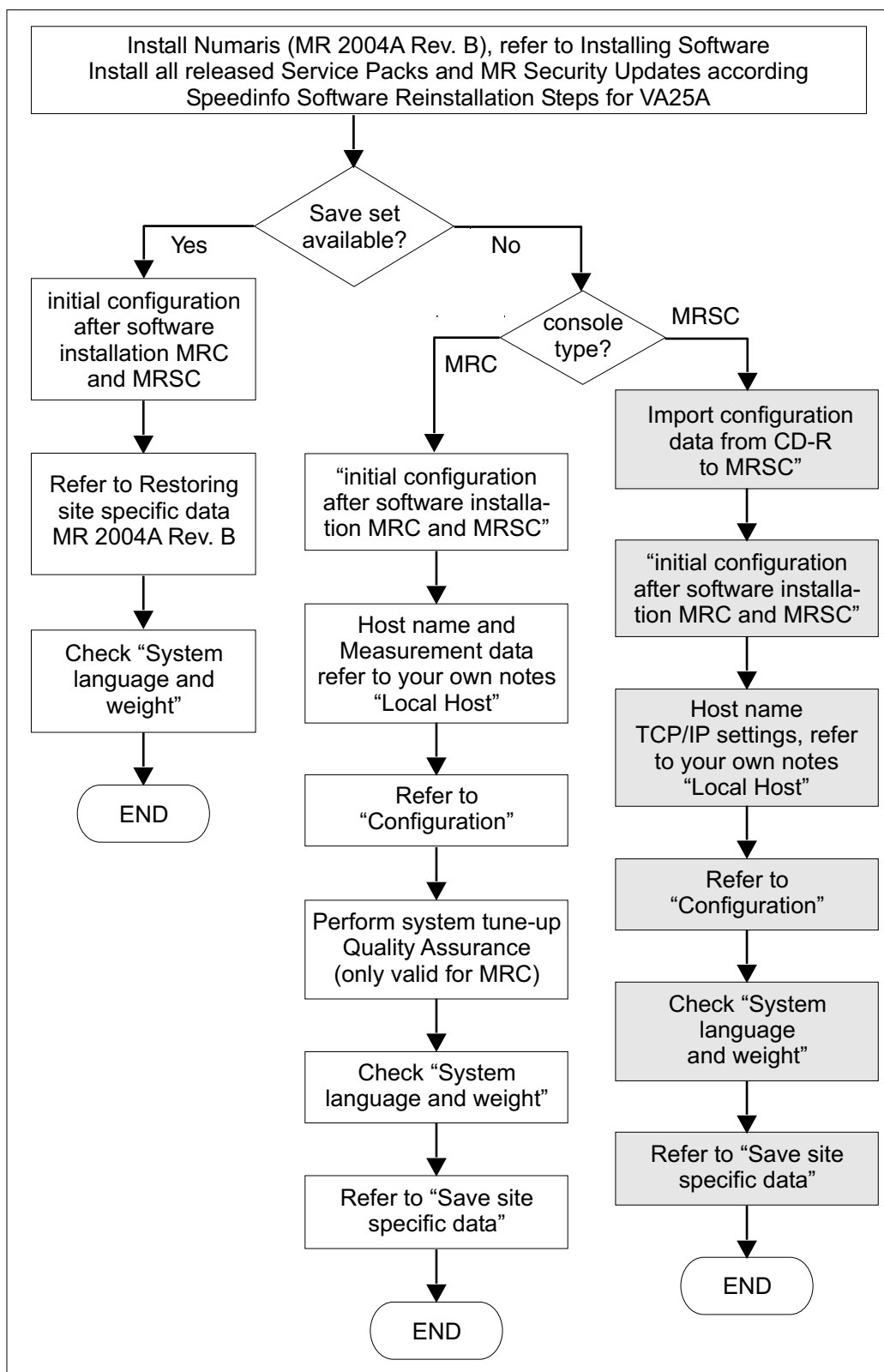
MPPS		
Input field		Your notes
Accession No.	Description	
Modality	Station	
Status	Type	
End Date	End Time	
Dose Dialog	Comments	

Worklist Results		
Input field		Your notes
Displayed HISRIS worklist result attributes		
Date of birth	Sex	
Accession Number	Admission ID	
Referring Physician	Performing Physician	
Requesting Physician	Ward	

3.1 Overview of installing VA25A MRC/MRC+MRSC

Recommendation: For paper versions, take out this page and use it for navigation.

Fig. 2: Overview of installing VA25A MRC/MRC + MRSC



3.2 General

The following chapter describes the software installation of VA25A on an empty system disk at the MRC and MRSC for MAGNETOM systems.



Make sure that you have all customer licenses as well as backups of site-specific data before reinstalling the software.

CV license filter - see (→ Requested data for cvlicense generation / Page 141).



The Shim files are part of the dynamic data.



In case of host replacement the local TCP IP settings are not restored.



If the system has suffered from HD crash refer to (→ Software / MR-000.840.11) TSG for troubleshooting information before you continue with the SW installation.

Table entries highlighted in gray are valid only for the MRSC installation:

This part of the table is always valid.
This part of the table is only valid for MRSC consoles

3.3 Installing Software

3.3.1 Prerequisites

- BIOS CD "MR059 BIOS CD VA25E with part no.: 100 49 286P9901"
- CD: 1/6 (International OS) 2/6 (XP Professional, syngo VD20N, CD 3/6 (NUMARIS/4 VA25A) is available.
- **(In case of MRC / MRSC host replacement including dongle. The "Dallas dongle" (Black) the application licenses are generated in relation to the Ethernet ID of the NIC. New application licenses must be generated in relation to the new Host)**



The application licenses are part of the save set (back-up)

Only if the back-up is missing you can use the original delivered Licence keys on floppy or CD disk.



The Service Licence is required for the configuration.

- Disconnect the LAN cable(s)



Consoles with the optional Pioneer 1.7 MOD, the power must be switched on for the MOD drive during the software installation.

- Enter the BIOS

R630
press F2
Enter password
Press F9 -Y
Insert BIOS CD "MR059 BIOS CD VA25E withpart no.: 100 49 286 P9901"
Boot Device Priority
move CD-ROM to the top (with the + key)
Press F10 - Y

- Wait until: "System Ready for Software Installation" is displayed.
- Remove the BIOS CD

3.3.2 Installing MR 2004A Rev. B MRC / MRSC

- Switch on the console to be installed (MRC or MRSC).
 - Insert the **(MR 2004A 1/6 International OS or CD5/6: Chinese OS or CD6/6: Japanese OS)** into the CD-ROM drive. **If the SW re installation fails e.g. due to a interrupt during the initial installation part of CD 1/6 refer to (→ Re installation of software fails / Page 120)**

Output	Input	Note
SELECTION WINDOW continue is highlighted	↵	Press any key to continue
	Automatic reboot ...	
SELECTION MENU Distribution type: Customer ... Create Partitions on ... CONTINUE INSTALLATION	Distribution Select CONTINUE INSTALLATION with the arrow key Press ENTER (↵)	
Confirm that "All DATA on selected harddisk) will be lost"	↵	Press ENTER (Installation time ~7 minutes)
INFORMATION WINDOW Please insert the next installation media labeled: syngo MR 2004A CD 2/6 Press any key to continue .. - Copying in progress	Please insert CD 2/6 Press ENTER (↵) (Installation time ~4 minutes)	
INFORMATION WINDOW Please insert the next installation media labeled: syngo MR 2004A CD 3/6 Press any key to continue .. - Copying in progress	Please insert CD 3/6 Press ENTER (↵) (Installation time ~4 minutes)	
The CD-ROM drive opens the carriage automatically for CD 3/6		Remove the CD from the CD-Rom drive

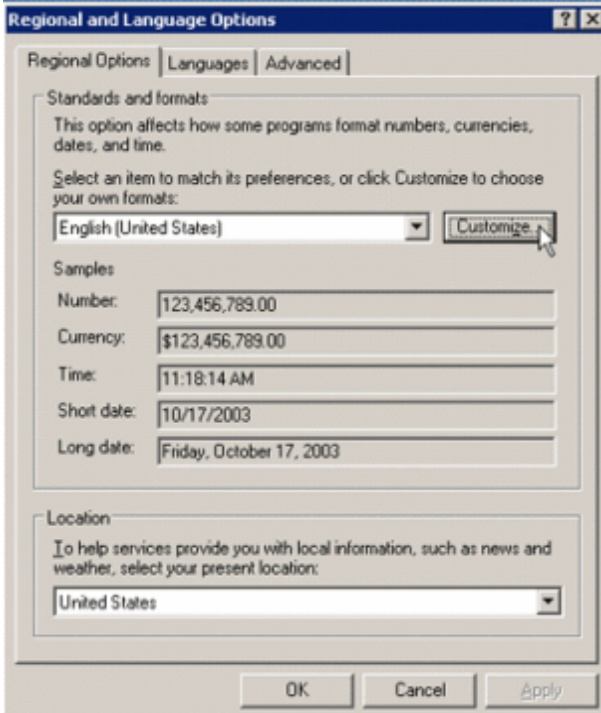
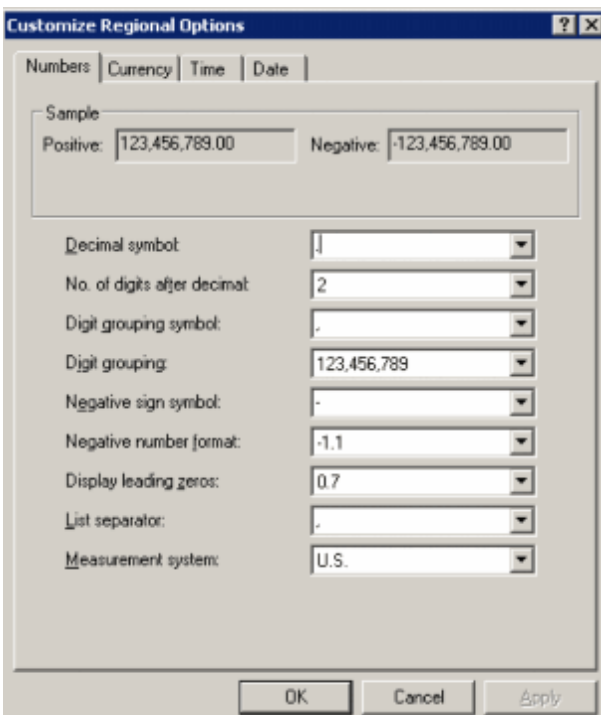
Output	Input	Note
Automatical reboot Windows XP Professional Setup Please wait while Setup copies files to your hard disk		Do not interfere the automatic installation / configuration. Windows XP Professional Setup will follow Several re-boots will follow Estimated time for the automatic configuration ~72 minutes
After installation is complete, the syngo MR 2004A Logon window is displayed.		
Install all released Service Packs and MR Security Updates according Speedinfo Software Reinstallation Steps for VA25A		

3.3.3 Keyboard settings and Digit grouping symbol for administrator

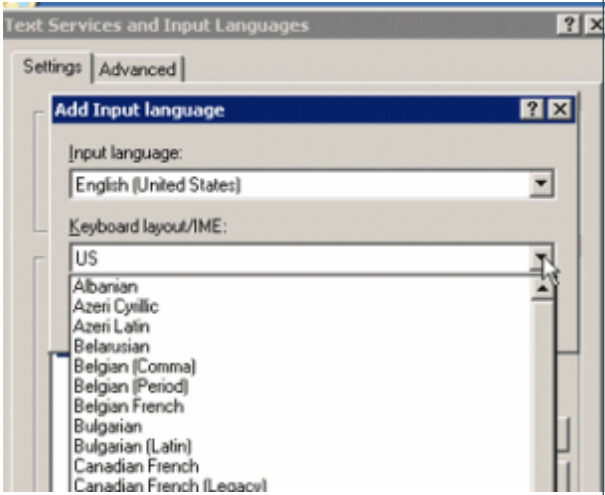
Following procedure is only to be performed for systems with **non-US** keyboard.

- Login as user **administrator**

Output	Input	Note
	Click START -> Settings -> Control Panel	Setting the local keyboard type
 Date, Time, Language, and Regional Options	Click Date, time, Language, and Regional options at the bottom right in the control panel	Only to be performed for systems with non-US keyboard
 Regional and Language Options	Click Regional and Language Options at the bottom right in the control panel	

Output	Input	Note
Fig. 3:  The 'Regional and Language Options' dialog box is shown with the 'Regional Options' tab selected. Under 'Standards and formats', the 'English (United States)' is selected in the dropdown menu, and the 'Customize...' button is highlighted with a mouse cursor. Below this, 'Samples' are listed for Number, Currency, Time, Short date, and Long date. The 'Location' section shows 'United States' selected in the dropdown menu. At the bottom are 'OK', 'Cancel', and 'Apply' buttons.	Click Customize	
Fig. 4: Decimal / Digit grouping symbols  The 'Customize Regional Options' dialog box is shown with the 'Numbers' tab selected. Under 'Sample', 'Positive' is '123,456,789.00' and 'Negative' is '-123,456,789.00'. Below this, various settings are listed with dropdown menus: 'Decimal symbol' (.), 'No. of digits after decimal' (2), 'Digit grouping symbol' (,), 'Digit grouping' (123,456,789), 'Negative sign symbol' (-), 'Negative number format' (-1.1), 'Display leading zeros' (0.7), 'List separator' (.), and 'Measurement system' (U.S.). At the bottom are 'OK', 'Cancel', and 'Apply' buttons.	Select the Numbers register Set "Decimal Symbol" to "." (= Full stop) and "Digit grouping symbol" to "," (= Comma). Click OK	

Output	Input	Note
Fig. 5: 	Click at the Languages, register; Click the Details Button Click Add	

Output	Input	Note
Fig. 6: 	Click pull-down menu "Keyboard layout/IME" Select input language related to your keyboard layout Click OK	
Fig. 7: 	Click pull-down menu "Select one of the installed input languages" Click Apply Click OK	The keyboard settings will be valid only for user administrator



If the system consists of MRC and MRSC check that both consoles are allowed to access each other through the ipsec filter.

The syngo-MR systems are protected with a service "ipsec".

Other network nodes, like MRSC or MRC and DICOM nodes, can only access the DICOM services / DB access that your console provides if they are included in the "Generic Transport Filter" of ipsec.

To list nodes that are allowed to access your console, click **Reports**; select **IP Security** below **SiteInfo**". To add a node manually into the IP security list, use following procedure:

At the "MRC", simultaneously press **Ctrl - Esc**. Select **Start, Advanced User**; enter the advanced user password; click **OK**

Simultaneously press **Ctrl - Esc**, click **Start, Run cmd** enter the command: **runas / user:administrator cmd** enter the password

In the newly started cmd tool you have administrator rights.

(replace the address "146.254.222.333" in the following **example** with the IP address of your remote node)

Input in the cmd window	Output
ipseccmd -f (0+146.254.222.333)	The command completed successfully

What to do next?		
Your Console type?		
MRC	MRSC with a valid save set	MRSC WITHOUT save set
Re-boot the console see (→ Reboot / Page 62) 		Keep login as administrator see (→ Export MRC configuration data to MRSC / Page 60) 

3.3.4 Export MRC configuration data to MRSC



If you have a valid save set for the MRSC skip this procedure. Continue with (→ Reboot / Page 62)

It makes sense only to perform the save at CD-R when you have completed the configuration of Customer UI: *syngo*, DICOM: Hardcopy devices and Network nodes at the MRC.

- Login as **administrator** at the **MRC**
- Insert an empty CD-R disk into the CD-R drive of the **MRC**

Output	Input	Note
	Click START Select RUN	Click START in the lower left corner Select RUN

Output	Input	Note
	cmd Click OK	
	mrscsave.bat ↵ Press Enter	

3.3.5 Import configuration data from CD-R to MRSC

Only syngo UI settings, DICOM configuration data and TCP/IP addresses from remote nodes are copied with this procedure.

Systems consisting of Main and Satellite consoles

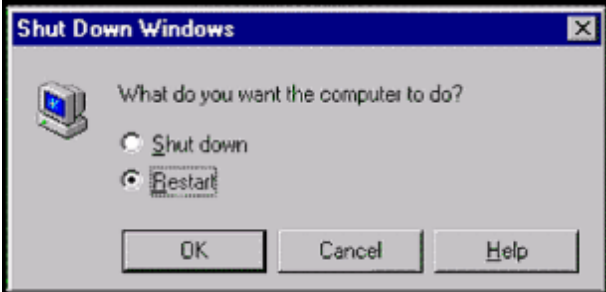
- You are logged in / login as **administrator** at the **MRSC**
- Insert CD-R disk into the CD-ROM drive of the **MRSC**, containing partial configuration data from the MRC (saved with the procedure "Export MRC configuration data to MRSC")

Output	Input	Note
	Click START Select RUN	Click START in the lower left corner Select RUN
	cmd Click OK	
	mrsc_default.bat ↵	



The MRSC will use the database of the MRC for export/import to remote DICOM nodes for shared DB configuration.

3.3.6 Reboot

Output	Input	Note
	Close all open windows	
Fig. 8: 	Click Start -> Shut Down Select Restart Click OK	Reboot the host

What to do next?
(→ Initial configuration after software installation MRC and MRSC / Page 63) 

3.4 Initial configuration after software installation MRC and MRSC

Leave the customer LAN disconnected



The "Ping" command does not work during the initial configuration since the LAN is disconnected.

- Login as user **meduser**



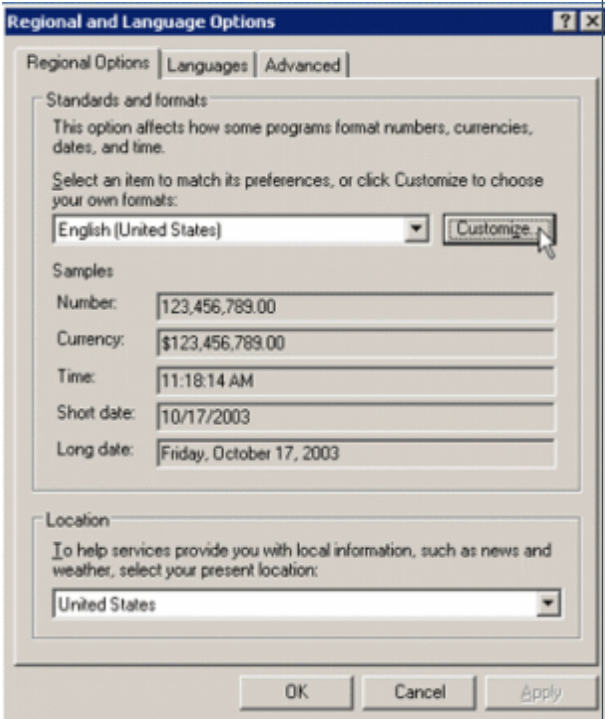
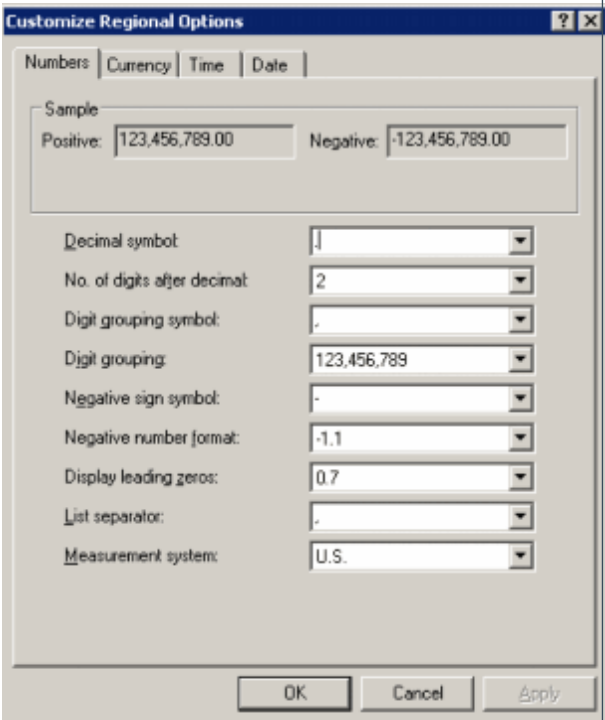
Before you enter into the displayed service platform, check/set keyboard settings.

3.4.1 Keyboard settings and Digit grouping symbol for meduser

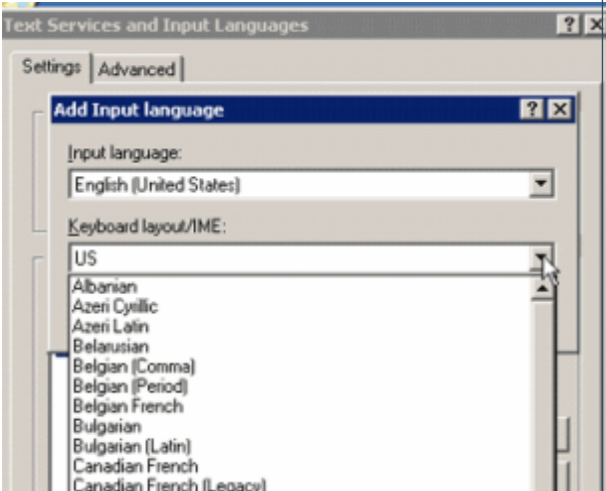
Following procedure is only to be performed for systems with **non-US** keyboard.

- Login as user **meduser**

Output	Input	Note
	Click START -> Settings -> Control Panel	Setting the local keyboard type
 Date, Time, Language, and Regional Options	Click Date, time, Language, and Regional options at the bottom right in the control panel	Only to be performed for systems with non-US keyboard
 Regional and Language Options	Click Regional and Language Options at the bottom right in the control panel	

Output	Input	Note
Fig. 9: 	Click Customize	
Fig. 10: <i>Decimal / Digit grouping symbols</i> 	Select the Numbers register Set "Decimal Symbol" to "." (= Full stop) and "Digit grouping symbol" to "," (= Comma). Click OK	

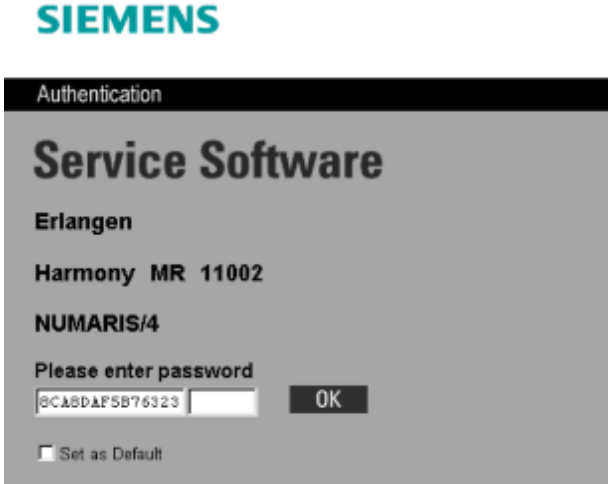
Output	Input	Note
Fig. 11: 	Click at the Registrar Languages, Click at Button Details Click Add	

Output	Input	Note
Fig. 12: 	Click pull-down menu "Keyboard layout/IME" Select input language related to your keyboard layout Click OK	
Fig. 13: 	Click pull-down menu "Select one of the installed input languages" Click Apply Click OK	The keyboard settings will be valid only for meduser

- 

Ignore message "Could not retrieve security data"
- 

Keep login during the complete initial configuration, DO NOT re-boot!

Output	Input
	Click the Local Service link to start configuration
Fig. 14:  The image shows the Siemens Service Software Authentication window. It has a title bar 'Authentication' and a main title 'Service Software'. Below the title, it says 'Erlangen' and 'Harmony MR 11002'. The user is identified as 'NUMARIS/4'. A prompt 'Please enter password' is followed by a text box containing '8C&BD&F5B76323' and an 'OK' button. At the bottom, there is a checkbox labeled 'Set as Default'.	Remember to click Set as default Enter the Service key (Level 7 key is required) Click OK

What to do next?		
Site specific data available? YES	Site specific data available? NO	
See "Restoring site specific data MR 2004A Rev. B on MR 2004A Rev. B" on Page 5 - 4 	Perform full configuration (Remember to import the customer licences) Refer to your own NOTES Chapter 2 "Getting started and configuration pages" (It's recommended to use the configuration pages for navigation purpose where you are in the configuration procedure) NOTE Information regarding known problems during the Initial configuration are described here separately for MRC and MRSC. 	
	Console type MRC: See (→ MRC / Page 68)	Console type MRSC: See (→ MRSC / Page 72)

3.5 MRC

For Help regarding configuration of the syngo MR system, please click **Help** at the top right of the Configuration platform refer to "Procedures MRC".

- After completed configuration click "Finish", re-boot the system as follows: Click "**Home**", confirm the re-boot message window with "**OK**", **please wait the system will need ~5 minutes before it will re-boot.**
- Reconnect the customer LAN connection



Check / install Hotfixes and Bugfixes

3.5.1 Activate BIOS protection

- Insert BIOS CD "MR059 BIOS CD VA25E with PN: 100 49 286 P9901"
- Reboot the system
- Wait until: "HIPAA installed" is displayed
- Remove the BIOS CD

3.6 IFN (Image Filter New) Installation Work Steps



The IFN is only available for newly delivered VA25A-based MR systems or Evolve Upgrades beginning with 2006 and up. In case of a software reinstallation, check if the following CD

"syngo MR 2004A,
Numaris/4 VA25A - IFN Rev. A,
CD 1/1,
Part no.: 10018855 K2219"

is available on the system. If not, the installation of the IFN can be skipped for the software reinstallation



For newly delivered VA25A-based MR systems or Evolve Upgrades beginning with 2006 and up, the IFN software replaces the image filter software from Contex Vision.

The filter function of the IFN (Image Filter New) is the same as that of the old one from Contex Vision.



The IFN can be installed only after the service pack 2 update.

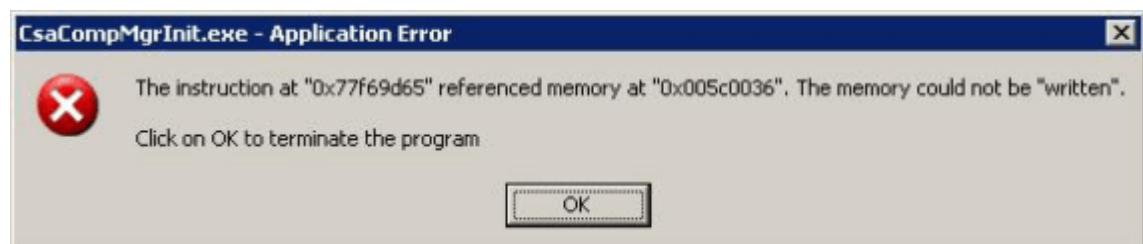


The IFN can be installed only before the service pack 3 update.

3.6.1 Hints

- On some systems, the following or a similar error message may be displayed occasionally:

Fig. 15:



3.6.2 Installing IFN on MRC and MRSC

ACTION	MRC	MRSC
- Check if SoftwareUpdate VA25A-SP2 is already installed. If not, please refer to the UI MR006/05/R or MR008/05/P. In the main syngo MR UI, you can verify if the SP2 is installed: - Help -> Info - Check if the "N4_VA25A_LATEST_20040724_P16" entry is visible under Siemens Software Packages or see logfile C:\MedCom\MriProduct\inst\SWHistory.txt	☐	☐
To install the IFN package, you have to be logged in as Administrator. If necessary, reboot and log in as Administrator (disable the autologin function in the local service UI).	☐	☐
Insert the CD "syngo MR 2004A, Numaris/4 VA25A - IFN Rev. A, CD 1/1, Part no.: 10018855 K2219" in the CD-ROM drive.	☐	☐
- Start Window Explorer: - Click START -> Programs -> Accessories -> Windows Explorer	☐	☐
Navigate to drive R:	☐	☐
Run program MR_Update.bat by double-clicking on it. NOTE: The program will first copy the software update package to the system. After copying is complete, installation of the package will start in a new CMD window (this will take approx. 20 - 30 min!) The original CMD window displays 'Press any key to continue ...' This can be done safely at any time. The CD will be ejected AUTOMATICALLY. Logfile entries such as: "Installation finished with ERRORLEVEL: 42" are indicating no errors!! The installation has finished when the following message is displayed in the CMD window entitled 'Software Update Installation - C:\sysmgmt\sd_store\MR_Update.bat INST': Finished software update at "actual date" "actual time", please reboot system!	☐	☐
Wait until please reboot system! is displayed.	☐	☐

ACTION	MRC	MRSC
Remove the CD "syngo MR 2004A, Numaris/4 VA25A - IFN Rev. A, CD 1/1, Part no.: 10018855 K2219" from the CD-ROM drive.	☐	☐
Open a cmd tool window.	☐	☐
Type in the following commands: SET TARGET=%Medhome%\MriProduct\ExamDbNls\en move /y %TARGET%\MrExamDb_EdbResen.dll %TARGET%\MrExamDb_EdbResen_old.dll copy /y %TARGET%\MrExamDb_EdbRes.dll %TARGET%\MrExamDb_EdbResen.dll	☐	☐
Reboot the complete system.	☐	☐

3.6.3 Verifying correct IFN installation

ACTION	MRC	MRSC
- In the main syngo MR UI, you can verify correct installation as follows: - Help -> Info - Check if the "N4_VA25A_LATEST_20040724_P25" entry is visible under Siemens Software Packages. If not, the IFN is not installed correctly.	☐	☐

3.6.4 CV License (only valid for systems delivered before 01/2006)

- If you have the option CV license filter, the low level key must be entered in the Service UI Utilities, Escape to OS with the command: **cvlicense** [license key]. For detailed information see (→ Requested data for cvlicense generation / Page 141).

What to do next?	
	Refer to the flowchart (→ Fig. 2 Page 52) after finishing the configuration

3.7 MRSC

For Help regarding configuration of the `syngo` MR system, please click **Help** at the top right of the Configuration platform refer to "Procedures MRSC"

3.7.1 General



An empty CD-R is used as medium for the procedure (→ Export MRC configuration data to MRSC / Page 60) and (→ Import configuration data from CD-R to MRSC / Page 61)
(not a floppy as described in the On-line Help texts)



SOFTWARE ERROR IMPORTANT!!!
When you perform Initial configuration
Database
Satellite Console
DON'T CLICK SAVE!
Refer to below workflow:

Workflow to avoid unwanted re-boot		
Local Host		
	System Type	
	Site Info	
	TCP/IP LAN	
	Users	
	MonitorType	
Database		
	Satellite Console	DON'T click SAVE just use > to navigate
Service		
	Licensing	Import your MRSC licenses
Navigate back to "Database"		

Workflow to avoid unwanted re-boot			
	Satellite Console	Perform configuration for your type of satellite console, save the system will re-boot automatically when you click save	
MRSC shares DB with MRC: See (→ MRSC Sharing Database with MRC / Page 73)		or	Local DB at MRSC: See (→ Completed configuration / Page 73)

3.7.2 MRSC Sharing Database with MRC

For Help regarding configuration of the syngo-MR system, please click Help at the top right of the Configuration platform refer to "Procedures MRSC"

3.7.3 Completed configuration

- After completed configuration click "Finish", re-boot the system as follows: Click "**Home**", confirm the re-boot message window with "**OK**", **please wait the system will need ~7 minutes before it will re-boot.**
- Reconnect the LAN connections



Check / install Hotfixes and Bugfixes

- Make sure that software options with additional installation media (CD's or DVD's) are installed

3.7.4 Activate BIOS protection

- Insert BIOS CD "MR059 BIOS CD VA25E with PN: 100 49 286 P9901"
- Reboot the system
- Wait until: "HIPAA installed" is displayed, remove BIOS CD from the CD

3.8 IFN (Image Filter New) Installation Work Steps



The IFN CD is only available for newly delivered VA25A-based MR systems beginning with 2006 and up. In case of a software reinstallation, check if the following CD

"syngo MR 2004A,
Numaris/4 VA25A - IFN Rev. A,
CD 1/1,
Part no.: 10018855 K2219"



is available on the system. If not, the installation of the IFN can be skipped for the software reinstallation



For newly delivered VA25A-based MR systems or Evolve upgrades beginning with 2006 and up, the IFN software replaces the image filter software from Context Vision. The filter function of the IFN (Image Filter New) is the same as that of the old one from Context Vision.



The IFN can be installed only after the service pack 2 update.



The IFN can be installed only before the service pack 3 update.

3.8.1 Hints

- On some systems, the following or a similar error message may be displayed occasionally:

Fig. 16:



3.8.2 Installing IFN on MRC and MRSC

ACTION	MRC	MRSC
- Check if SoftwareUpdate VA25A-SP2 is already installed. If not, please refer to the UI MR06-05-R or MR08-05-P. In the main syngo MR UI, you can verify if the SP2 is installed: - Help -> Info - Check if the "N4_VA25A_LATEST_20040724_P16" entry is visible under Siemens Software Packages or see logfile C:\MedCom\MriProduct\inst\SWHistory.txt	☐	☐
To install the IFN package, you have to be logged in as Administrator. If necessary, reboot and log in as Administrator (disable the autologin function in the local service UI).	☐	☐
Insert the CD "syngo MR 2004A, Numaris/4 VA25A - IFN Rev. A, CD 1/1, Part no.: 10018855 K2219" in the CD-ROM drive.	☐	☐
- Start Window Explorer: - Click START -> Programs -> Accessories -> Windows Explorer	☐	☐
Navigate to drive R:	☐	☐
Run program MR_Update.bat by double-clicking on it. NOTE: The program will first copy the software update package to the system. After copying is complete, installation of the package will start in a new CMD window (this will take approx. 20 - 30 min!) The original CMD window displays 'Press any key to continue ...' This can be done safely at any time. The CD will be ejected AUTOMATICALLY. Logfile entries such as: "Installation finished with ERRORLEVEL: 42" are indicating no errors!! The installation has finished when the following message is displayed in the CMD window entitled 'Software Update Installation - C:\sysmgmt\sd_store\MR_Update.bat INST': Finished software update at "actual date" "actual time", please reboot system!	☐	☐
Wait until please reboot system! is displayed.	☐	☐

ACTION	MRC	MRSC
Remove the CD "syngo MR 2004A, Numaris/4 VA25A - IFN Rev. A, CD 1/1, Part no.: 10018855 K2219" from the CD-ROM drive.	☐	☐
Open a cmd tool window.	☐	☐
Type in the following commands: SET TARGET=%Medhome%\MriProduct\ExamDbNls\en move /y %TARGET%\MrExamDb_EdbResen.dll %TARGET%\MrExamDb_EdbResen_old.dll copy /y %TARGET%\MrExamDb_EdbRes.dll %TARGET%\MrExamDb_EdbResen.dll	☐	☐
Reboot the complete system.	☐	☐

3.8.3 Verifying correct IFN installation

ACTION	MRC	MRSC
- In the main syngo MR UI, you can verify correct installation as follows: - Help -> Info - Check if the "N4_VA25A_LATEST_20040724_P25" entry is visible under Siemens Software Packages. If not, the IFN is not installed correctly.	☐	☐

3.8.4 CV License (only valid for systems delivered before 01/2006)

- If you have the option CV license filter, the low level key must be entered in the Service UI Utilities, Escape to OS with the command: **cvlicense** [license key]. For detailed information see (→ Requested data for cvlicense generation / Page 141).

What to do next?	
	Remember to create a backup if the configuration has changed (refer to (→ Save site specific data / Page 80)

4.1 Overview installing Imager SW VA25

4.1.1 General

This chapter describes how to install the Imager software.



The MAGNETOM syngo-MR Imager is delivered with factory-installed software. The following Installation procedure is necessary only in case of replacement of the host.



Before performing one of the following steps, refer to (→ Host / Imager / MR-000.840.12) for troubleshooting information.

4.2 Installing Software at the Imager "R630"

4.2.1 Prerequisites

- Customer-CD 4/6 MR 2004A Rev. B
- The measurement network is connected

4.2.1.1 Check/Import the Imager license to the MRC

- Login as administrator (MRC)

Use Windows file Explorer to copy/paste the license file from the floppy disk to the following target on the MRC: **C:\MedCom\tftp**

- Rename the license file for the Imager to **license.dat**

4.2.2 Installation procedure

- Insert the CD "Customer-CD 4/6." in the CD-ROM drive
(Windows XP, syngo, and syngo MR is included at CD 4/6)
- Switch off / on the power to force a re-boot from the CD-ROM drive
- Please wait ~10 minutes, the CD-ROM drive opens the carriage automatically, remove the CD.
- The silent installation will need ~45 minutes additionally

Following table is only provided as additional information for sites with a monitor and keyboard connected at the imager:

Output	Note	
A pseudo-graphic surface appears with the warning that all data will be lost after pressing a key.		All data will be lost at the system disk of the Imager. (Installation time ~10 minutes)
The Imager will reboot		
The CD-ROM drive opens the carriage automatically for CD 4/6		Remove the CD from the CD-ROM drive
Automatic reboot Windows XP Professional Setup Please wait while Setup copies files to your hard disk	Do not interfere the automatic installation / configuration. Windows XP Professional Setup will follow Several re-boots will follow Estimated time for the automatic configuration ~ 45 minutes	



The Imager license must be available at the MRC

- Boot the MRC (Host).
- Start the MR application.
- Switch on power to the Imager



The Imager will reboot twice following the installation.

The initialization of each raw-data disk takes approximately 5 minutes!

(No activity will be seen in the task manager)

(To check the measurement network connectivity to the Imager enter the command **ping -t mrir** in a cmd tool at the MRC)

- The imager is ready when: on each image drive F:\ (at High-End extension additionally G:\, H:\, I:\) there is a file with the size of 4 GB placed with the name "rawex-change.dat".

The imager is ready to use, when this file exists on all the image drives with the above-mentioned size.

(Additional files at F: "MattauchBuffer1.dat" and G: "MattauchBuffer2.dat") (You can use the Remote Desktop Connection to access the Imager: Start, Programs, Accessories, Communications, Remote Desktop Connection. The computer name of the Imager is **mrir**, User: **MRSystem**).

What to do next?

Install all released Service Packs according Speedinfo Software Reinstallation Steps for VA25A

Make sure that the system is producing images.

5.1 Save site specific data

Save the site-specific data each time you make changes to the site-specific configuration, tune-up, or sequences and protocols. In case of a hard disk crash, you will be able to restore a valid set of site-specific data.

Perform a backup of the following packages:

5.1.1 Numaris4

The Numaris4 package contains the complete system tune-up data, e.g. shim files, coil tune-up, measurement system-specific data

5.1.2 CustomerProtocols

The CustomerProtocols package contains the customer protocols.

5.1.3 SW-Settings

The SW-Settings02 package includes configuration data regarding e.g. remote nodes, DICOM nodes, Licensing, Customer UI...

- Table entries highlighted in gray are valid only in certain cases. Refer to the "Note".

This part of the table is always valid.
This part of the table is only valid for special cases. Refer to the related "Note".



Backups to local disks can be performed only at the customer site, not via remote logon.

5.1.4 Prerequisites

You need a empty recordable compact disk medium.



Only empty CD-R (recordable compact disk) media may be used to save site-specific data.



If there is a HC Device "DSS" in the backup data set, the restore of the camera settings will fail.

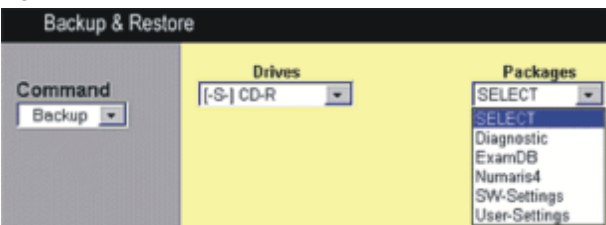
5.1.5 Settings that are not part of the backup/restore

The following settings are not part of the site specific data. Make sure to note the settings in your configuration pages.

- Local Time zone
- **Local host, Users, "Enable Autologin"**
- Basic Worklist SCP Query Properties, Query Waiting Time
- **Printer installation, printer configuration**
- Static routes (See: (→ Check / add route to SIEMENS Remote Services router / Page 100))
- System language and weight
- **Settings for "Remote Service Access Control** (e.g., "Permanent Limited Access")

5.1.6 Backup (Save site specific data) on CD

- Login as **administrator** ↵.
- Start the service platform, select, **Backup & Restore**.

Output	Input	Note
Fig. 17: 	Select Backup in the Command pull-down menu	Select Backup
Fig. 18: 	Select the backup device in the Drives pull-down menu, e.g. CD-R	The CD-Burner has the drive letter S :
CustomerProtocols	Select CustomerProtocols in the Packages pull-down menu Click Go	Estimated backup time is ~ 1-12 minutes
SW-Settings02	Select SW-Settings02 in the Packages pull-down menu Click Go	Estimated backup time is ~ 1 minutes
Wait for message: "End Backup"		

Output	Input	Note
Fig. 19: Numaris 4 	Select Numaris4 in the Packages pull-down menu Click Go	Estimated backup time is ~ 3-15 minutes
	Wait for message: "End Backup"	
SW-Settings02	Select MriCVLicense-File in the Packages pull-down menu Click Go	Estimated backup time is ~5 seconds
	Wait for message: "End Backup"	

5.1.7 Backup (Save site specific data) on USB device



Please check the external media (e.g. USB-disk) first with the virus scanner of the service PC. Be sure that the external media used is free of viruses.

- Login as **administrator** ↵.
- Start the service platform, select, **Backup & Restore**.
- Start the "Internet Explorer"
- Connect the USB device.

Output	Input	Note
	Select the USB device in the Drives pull-down menu, e.g. F	E.g. drive letter F:
CustomerProtocols	Select CustomerProtocols in the Packages pull-down menu Click Go	Estimated backup time is ~ 1-8 minutes
SW-Settings02	Select SW-Settings02 in the Packages pull-down menu Click Go	Estimated backup time is ~ 1 minutes
	Wait for message: "End Backup"	
SW-Settings02	Select Numaris4 in the Packages pull-down menu Click Go	Estimated backup time is ~ 2-10 minutes
	Wait for message: "End Backup"	

Output	Input	Note
SW-Settings02	Select MriCVLicense-File in the Packages pull-down menu Click Go	Estimated backup time is ~5 seconds
	Wait for message: "End Backup"	

5.1.8 Additional packages

Output	Input	Note
The following packages have to be backed up in certain cases only:		
Security-Settings Only valid for systems with HIPAA	Select Security-Settings in the Packages pull-down menu Click Go	Only to be performed at the MRC (and stand-alone MRSC's with a Local DB) . Security and user / group settings of HIPAA configuration
Diagnostic Performed only in special service cases!	Select Diagnostic in the Packages pull-down menu Click Go	Files generated with MrSaveLog.bat Estimated back-up time is less than 1 minute
ExamDB Only valid for application specialists for special cases!	Select ExamDB in the Packages pull-down menu Click Go	Only to be performed at the MRC Estimated backup time 1-25 min, depending on the number of protocols

5.2 Restoring site specific data MR 2004A Rev. B on MR 2004A Rev. B

Restoring the site-specific data may be helpful in many situations. Please read the following information before beginning work.

- Troubleshooting: In case of corrupted site-specific data, we recommend saving the actual state of the site-specific data to another medium before restoring an old save set. This will save time if the problem is resolved by reinstalling the software. The actual tune-up data may save a great deal of work, especially if the customer recently optimized his protocols.
- Disk crash: After installing the software, a set of site-specific data may prove to be helpful.



Backups to local disks can be performed only at the customer site, not via remote login.

5.2.1 Prerequisites.

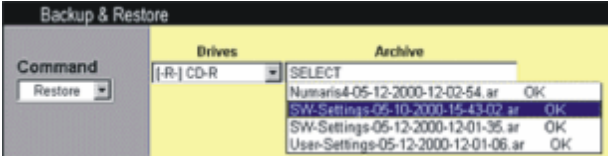


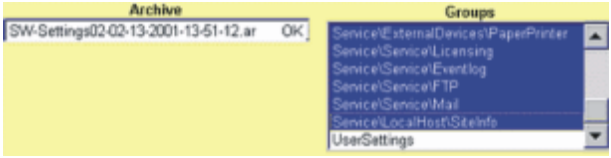
Ensure that the save set contains data valid for the current software version

5.2.2 Restore site-specific from CD.

(In case of initial configuration, a security warning is displayed. Confirm with "Yes" and use **Alt-Tab** to toggle to the Local Service Platform)

- Login as user: **administrator**
- Insert Backup CD-R, Select, **Backup & Restore** in the service platform

Output	Input	Note
Fig. 20: 	Select Restore in the pull-down menu Command Select device, e.g. CD-R the pull-down menu Drives	The CD-ROM will get the drive letter R :
Fig. 21: 	Select SW-Settings02 in the Archive pull-down menu	

Output	Input	Note
Fig. 22: 	Select all groups except "User-Settings and without SecConfig)" in the Groups pull-down menu	Due to a software bug, the User-Settings group must be restored separately after a system reboot
General use of mouse and key combinations To select all groups at once, click the group at the top, press and hold the Shift key, and click the group at the bottom. All groups are highlighted in blue. (If you press and hold the Ctrl key, any combination of groups can be selected by clicking the mouse)		
	Click Go	

- Wait until the message "***** End Restore:" is displayed

Output	Input	Note
	Select MriCVLicenseFile in the Archive pull-down menu	Low level license of Context Vision filter
	Click Go	

- Wait until the message "***** End Restore:" is displayed

Output	Input	Note
	Select Numaris4 in the Archive pull-down menu	Ignore error message "License file ...Error"
	Select all Groups	Confirm the popup window: 'Errors in RestoreHook' with OK
	Click Go	

- Wait until the message "***** End Restore:" is displayed
- Wait until the message "***** End Restore:" is displayed

Output	Input	Note
	Select CustomerProtocols in the Archive pull-down menu	Only restore the package "CustomerProtocols" if it's more recent then the package "Numaris4" (The package Numaris4 contains the customer protocols as-well.)
	Select all Groups	
	Click Go	

- Wait until the message "***** End Restore:" is displayed

- Click "**Home**" wait for the report generation to finish. Confirm pop-up with "**Cancel**"

5.2.3 Restore site-specific from external USB device



Please check the external media (e.g. USB-disk) first with the virus scanner of the service PC. Be sure that the external media used is free of viruses.

Restoring the site-specific data may be helpful in many situations. Please read the following information before beginning with the work.

In case of corrupted site-specific data, we recommend saving the actual status of the site-specific data to another medium before restoring an old save set. This will save time if the problem is resolved by reinstalling the software. The actual tune-up data may save a great deal of work, especially if the customer optimized his protocols recently.

- Login as **administrator** ↵.
- Start the "Internet Explorer"
- Login to the service platform, select **Backup & Restore**.
- Connect the USB device.
- Select **Restore** in the "Command" drop-down menu.
- Select the source USB-drive (drive letter e.g. **F:**) in the "Drives" menu.
- Select your packages to be restored xxx_xxx_xxxx_xxxx_ok.ar in the "Archive" menu.
- In the menu "Groups", select the necessary packages to be restored.
- Click **Go** to start the "Backup", check for errors in the log window after restore has been completed.



To ensure the Administrator account password is correct following restoration of the dynamic data, it is recommended to set the default Administrator password in the Service GUI > **Configuration** > **Local Host** > **Users**

5.2.4 Settings that must be manually activated after restore

Following configuration data needs manual actions to be activated after a restore has been performed.

- Local Time zone
- In case of host replacement the local TCP IP settings are not restored.
- **Local host, Users, "Enable Autologin"**
- Basic Worklist SCP Query Properties, Query Waiting Time
- **Printer installation, printer configuration**
- Static routes (See: (→ Check / add route to SIEMENS Remote Services router / Page 100))
- System language and weight

- **Settings for "Remote Service Access Control** (e.g., "Permanent Limited Access")



In case the error message "C:\MedCom\MriSitedata\ContextVision\.siemensIRU-PCNT-init-> access denied" is displayed during restore of N4 data, re-enter the cvlicense as described in (→ Context Vision / Page 141).

- Activate the measurement system configuration data as follows: Start the local Service platform, select **Configuration**, and click **"Save"** in each platform listed under **"measurement"** (the changes will be effective after the next reboot).
- Reboot the system from the **start** menu, select **Start -> Shut Down** -> select Re-start and press **"OK"**.
(in case of SW reinstallation, the reboot takes ~ 20 additional min.)
- Login as **meduser**
- Select **"SW-Settings02"** in the Archive pull-down menu. Select **"User-Settings"** in the Groups menu.
- Click **"Go"** to start the restore. Read the messages in the log window.
- Wait until the message "***** End Restore:" is displayed.
- Click **"Home"** wait for the report generation to finish.
- Confirm pop-up window with **"Cancel"** and then reboot system manually after closing the Service Platform - Reboot the system from the **start** menu, select **Start -> Shut Down** -> select Shutdown and press **"OK"**.
- **Re-boot the system** (Connect the customer LAN if disconnected)
- After the reboot, perform a spot check regarding image quality. Also, verify that all customer options are listed in the system imprint.



The HIPAA configuration is only to be restored at the MRC (and standalone MRSC's with a Local DB).

Output	Input	Note
Before restore of security settings a reboot of the Host is necessary.		

Output	Input	Note
		Login as meduser
Security-Settings Only valid for systems with HIPAA	Select Security-Settings in the Packages pull-down menu Click Go	Only to be performed at the MRC (and stand-alone MRSC's with a Local DB) . MRSC that shares the DB with MRC also shares the security settings with the MRC

6.1 General

For systems that include a satellite console, **configure the main console first.**



You need your notes from the configuration pages.



The configuration pages contain important information regarding the default configuration/settings for syngo MR systems refer to (→ Configuration pages for MRC / Page 20), for MRSC refer to (→ Configuration pages MRSC / Page 40).



MRC/MRSC are used as abbreviations for the main/satellite console, mrc/mrsc are used as host names.



Enabling of security settings (HIPAA) is to be performed by qualified personnel with HIPAA training.



The use of non-released combinations of hardcopy cameras and BU products is prohibited and does not meet the requirements defined. For detailed information regarding the CAMERA configuration, refer to technical documentation "Camera Info Numaris 4 VA25A"



The "SIEMENS Remote Services TM" will only support router connections for NEWLY delivered systems. See <http://cs.med.siemens.de/> "For Service, Siemens Remote Services" for the most recent information regarding router configuration.



If you have the option "Image Filter Software"(Context vision) an additional license must be enabled.

Refer to (→ Context Vision / Page 141).



For known bug, please refer to "Known bugs and hints for Service"



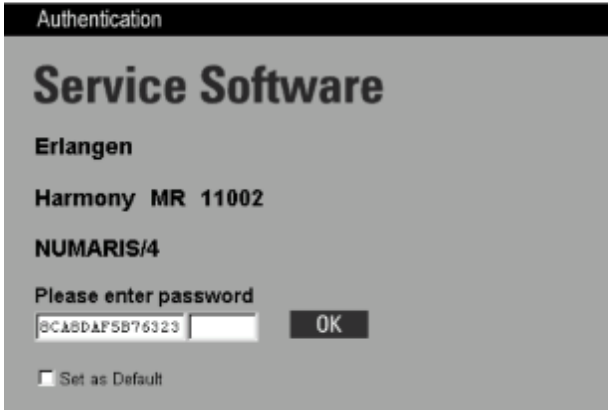
It's recommended **not** to configure WINS or DNS for **syngo** MR systems. (Netstat -a returns putpnt ~1h when WINS and DNS are configured)

6.2 Prerequisites

- The configuration pages contain the information/notes required for your system.
- The VA25A software is installed.

6.3 Access Service Software

- Simultaneously press the **Ctrl - Alt - Del** keys to open the login window.
- Login as user **meduser**.
- Click **Options** and select **LocalService**.

Output	Input
Security warning	Click Yes
Fig. 23: 	Remember to click Set as Default Enter the variable key Click OK
	Click Next



For detailed information regarding the configuration of syngo MR systems, refer to "Help" in the Service Software Platform.

6.3.1 Changes in the configuration

- If you have changed the configuration, re-boot the system as follows: Click "Home", confirm the re-boot message window with "OK"

What to do next?	
	Remember to create a backup if the configuration has changed (refer to (→ Save site specific data / Page 80)).

6.4 Parallel Printer Installation



System with two consoles (MRC and MRSC), local printer connected via the parallel port is only supported at the MRSC.

- Login as user Administrator

6.4.1 Installing local parallel printer

- Select Start, **Settings, Control Panel**.
- Click **Printers and Other Hardware**
- click **Printers and Faxes**
- Click the **Add Printer** icon.
- Click **Next** in the Add printer Wizard "Welcome to Add Printer Wizard"
- Select radio button **Local printer attached to this computer**, click **Next**
- Ignore that the Add Printer Wizard "Searching... " terminates with the message "The wizard was unable to detect any Plug and Play printers..."
- **Select LPT1 ... in the pull-down menu "Use the following port:", click Next**
- Select your printer in the Add printer Wizard "Install Printer Software", click **Next**
- Select radio button **Keep existing driver (recommended)**, click **Next**
- Name the printer in the Add printer Wizard "Name Your Printer", click **Next**
- Give a logical Name for the printer in the Add printer Wizard "Printer Sharing", click **Next**
- Select radio button **Yes** in the Add printer Wizard "Print Test Page", click **Next**
- Click the **Finish** in the printer Wizard "Completing the Add Printer Wizard".
- Click **OK** in the pop-up "A test page is now being sent to the printer ..."

6.4.2 Enabling MRC to access printer share of MRSC



If the actual console is MRSC, you must grant access to the MRC in the IP security filter.

- Keep login as user Administrator (at the MRSC)
- Open a cmd window
- enter the command **ipscmd -f (0+146.254.111.222)** exchange 146.254.111.222 in this example with the IP address of your MRC

6.5 Network Printer

The following procedure must be performed at each MR console that will print to the shared network printer.

- Login as user Administrator

6.5.1 Installing network printer

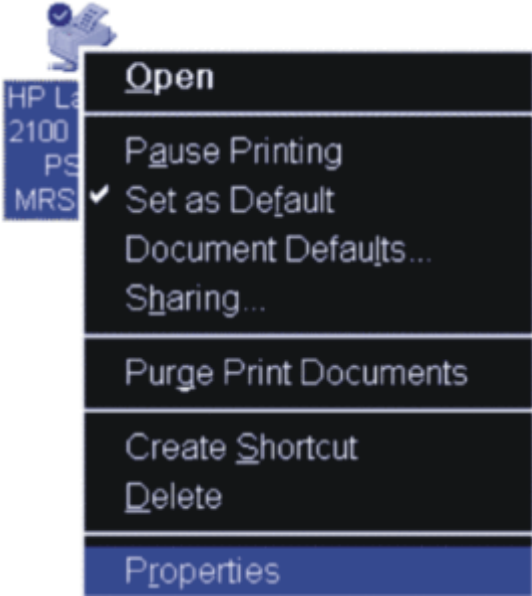
- Select Start, **Settings, Control Panel**.
- Click **Printers and Other Hardware**
- click **Printers and Faxes**
- Click the **Add Printer** icon.
- Click **Next** in the Add printer Wizard "Welcome to Add Printer Wizard"
- Select radio button **A network printer, or a printer attached to another computer**, click **Next**
- **Select radio button Connect to this printer (or to browse for a printer ... , enter the IP address of the MRSC e.g. \\146.254.111.221 in field "Name:" (In the Add printer Wizard "Specify a Printer"), click Next**
- Click **Yes** in the pop-up "Connect to Printer" You are about to connect to a printer on xx.xx.xx.xx ...
- Click the **Finish** in the printer Wizard "Completing the Add Printer Wizard".

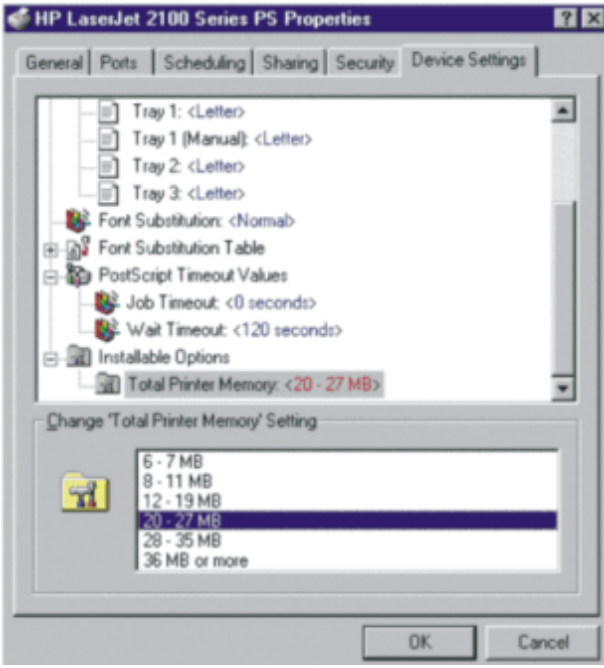
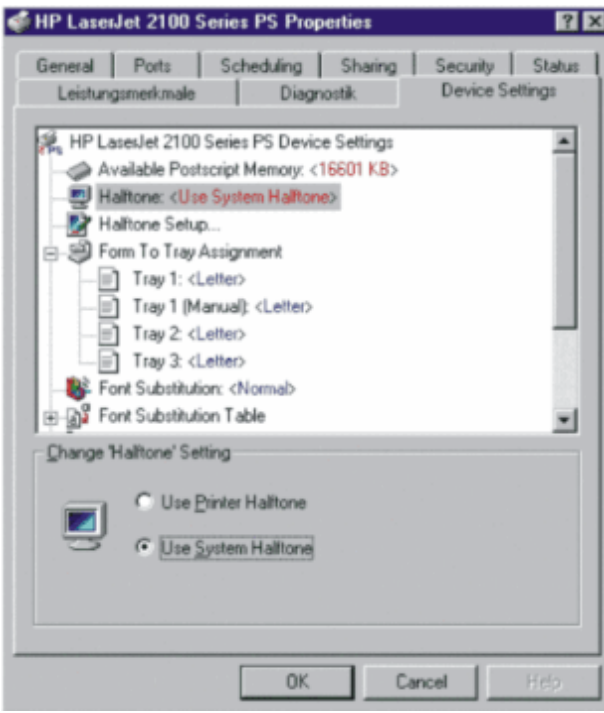
6.6 General Printer Properties

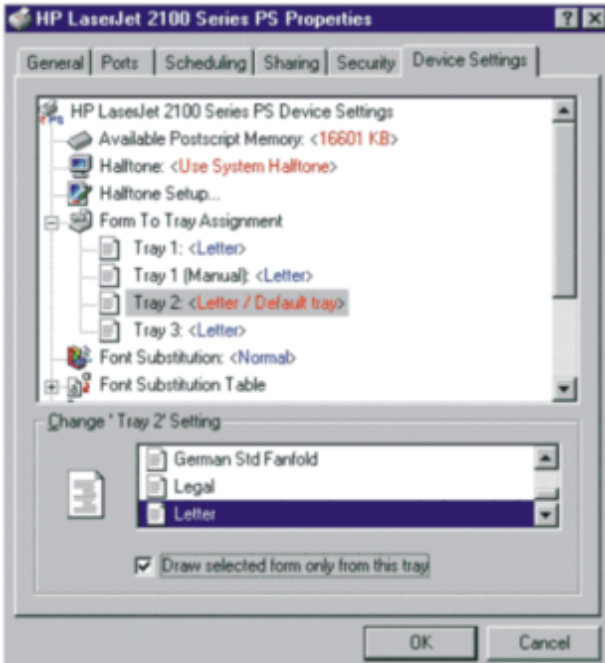
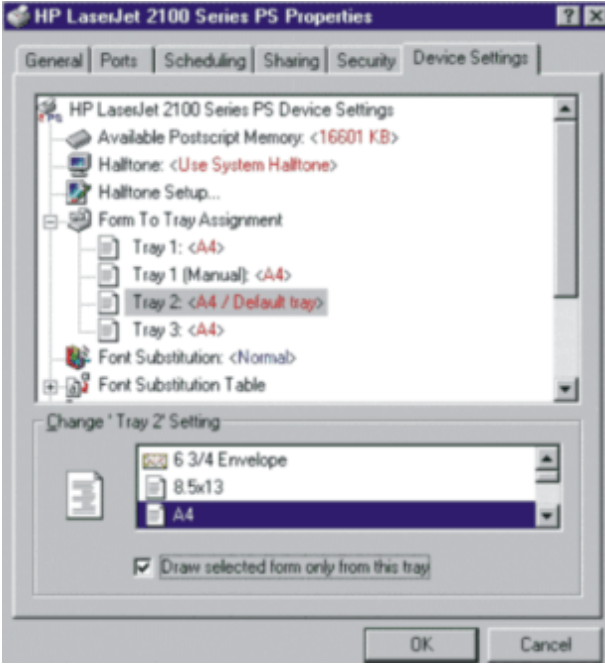
Each MR console that will be used to print has to be configured as follows to optimize print quality.

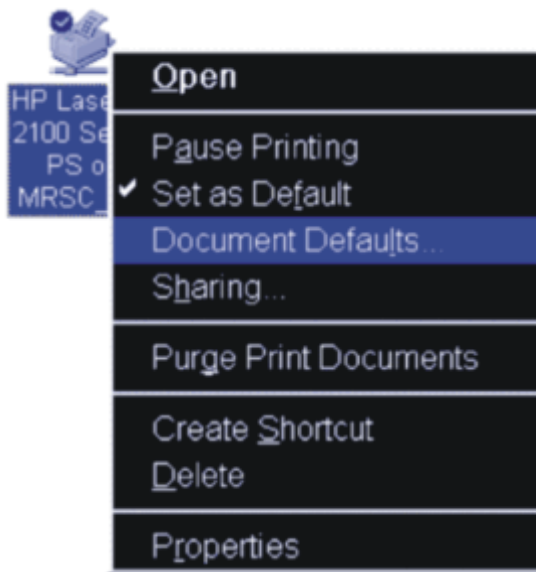
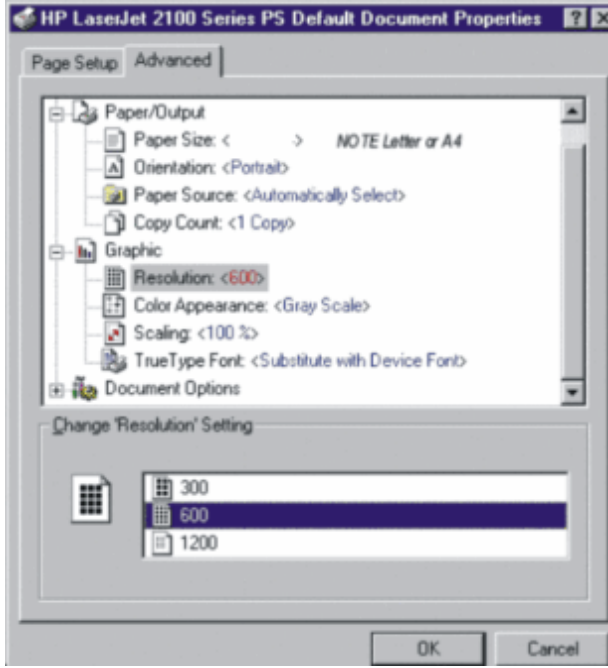


The printer should have at least 24MB RAM.

Output	Input	Note
	Click Start and select Settings, Printers	Click START in the lower left corner
Fig. 24: 	Right -click the HP printer icon and select Properties	

Output	Input	Note
Fig. 25: 	Click the Device Settings tab Navigate to Total Printer Memory Adjust the "Total Printer Memory Setting"	A minimum of 24 MB are recommended
Fig. 26: 	Check that half tone is set, select "Use System Halftone" radio button	

Output	Input	Note
Fig. 27: Paper size: US Letter 	Navigate to "Form To Tray Assignment" Click Tray2: (Select used paper size for the custom-er) Click the Draw selected form only from this tray check box	Only valid for printers using the "Letter" paper size
Fig. 28: Paper size: A4 	Navigate to "Form To Tray Assignment" Set each tray to A4 Click Tray2: Click the check-box Draw selected form only from this tray check box	Only valid for printers using the "A4" paper size
	Click OK	

Output	Input	Note
Fig. 29: Document Defaults 	Right-click the HP printer icon and select Document Defaults	
Fig. 30: Paper size & Resolution 	Click the Advanced tab Letter or A4 For A4 paper size, select A4 For "Letter", no changes are necessary Navigate to "Graphic" Select "Resolution" Set the resolution to 600 dpi	
		Click OK
Reboot the system		

- Remain logged in or log in as user "administrator"
- Log in as user "meduser".
- Start the local service platform, select **paper printer** under "external devices" and set the paper printer media formats available for your printer.

- Print a test image from the application to verify the print quality.

7.1 General

Only IE 6 has been tested and released for remote service as service browser at syn-go MR systems running with NUMARIS MR 2004A Rev. B.

The project manager is responsible to plan for a router connection to the Siemens Remote Diagnostic server. Regarding questions to your pre configured router from Siemens please contact (→ mailto:srs_final@med.siemens.de, mailto:srs_final@med.siemens.de) or your remote diagnostic responsible.

7.2 Prerequisites

- Router is pre-configured and known as route to SRS by the MRC (MRSC if available)
- System is configured at the SRS server

7.2.1 System Management

System management is a new remote service for syngo MR systems, it provides following functionality: proactive event reporting and statistical data transfer. In the future it will be used for SW distribution e.g hotfix / buxfix and anti virus patterns.

7.2.2 Pro Active Service

With MR 2004A Rev. B we provide Autoreports for proactive event triggering, that can be distributed on demand.

7.2.3 Quality Assurance DB (QDB)

The QA data will be transferred to the server at regular intervals for statistic evaluation and trend analysis, the data is collected into the MR-quality data-base.

Helium level and magnet shield temperature data is transferred to the server at regular intervals. This function offers a quick overview of the current He levels and magnet shield temperatures for the reporting sites.

Under following URL you can see data that has been sent by your system: <http://td-quality.med.siemens.de/>

7.2.4 Event DB (EDB)

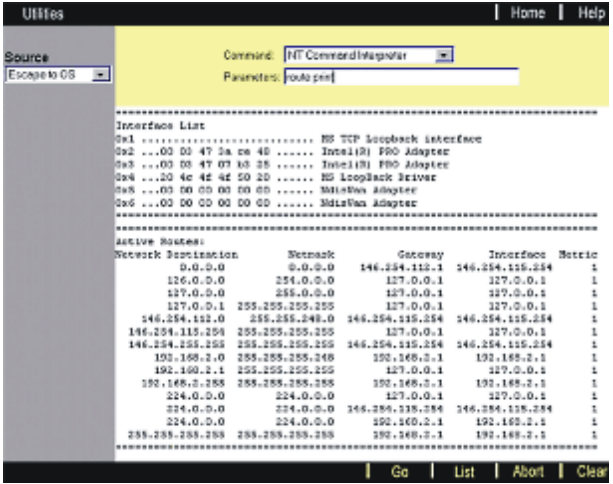
The EDB provides reports regarding events that are stored in the system event log.

7.2.5 System Utilization

Utilization Management is an new service offering and for leagal aspects the customer has to sign a contract in advance. Please involve your service organization before you activate System Utilization.

7.2.6 Check / add route to SIEMENS Remote Services router

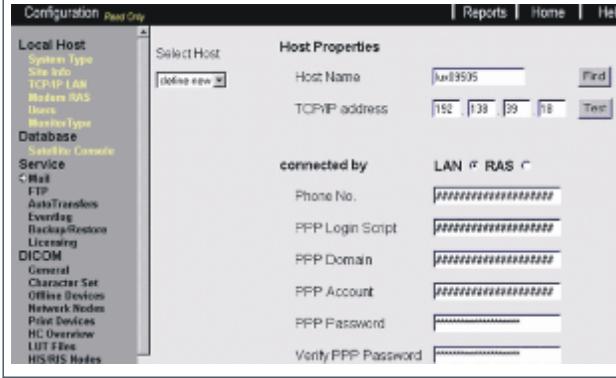
As the most optimal solution, the customer network administrator has configured the route from the default gateway to the SIEMENS SRS server for your continent. In the case of a VPN connection being used for SRS, the customer network administrator must configure the route from the default gateway to the SIEMENS SRS server.

	Output	Input
Fig. 31: Route print 		Start Local Service Select Utilities Select Escape to OS in the pull-down menu "Source" Enter the command in the Parameters field route print ↵ Check if the router for rdiag is listed In case the route is not listed to your SRS subnet you must add a static route.
Fig. 32: Route add 		The below example means use the customer router / gateway with IP: xx.xx.xx.xx for the SRS network as destination. E.g. (LUX09505 194.138.39.18) 194.138.39.0 route to the whole SRS subnet, and -p makes the entry persistent)
(Replace the xx.xx.xx.xx with your customer router / gateway IP address !!!)		
Europe, Africa: North-/ South America: Asia:	route add -p 194.138.39.0 mask 255.255.255.0 xx.xx.xx.xx ↵ route add -p 129.73.116.0 mask 255.255.255.0 xx.xx.xx.xx ↵ route add -p 194.138.243.0 mask 255.255.255.0 xx.xx.xx.xx ↵	

Check to ensure that configured SIEMENS Remote Server is added to the IP security list, click Reports; select IP Security below "SiteInfo". If the System Management server is missing, refer to (→ Adding host to IP security list / Page 139).

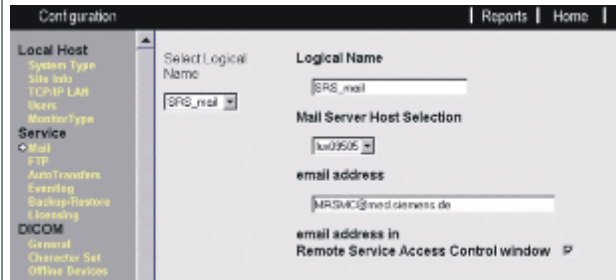
Only in case you have a wrong route, use the following command to delete: route delete 1x.x.x.x ↵	(To delete a wrong route from the routing table: enter IP address of the routing network)
--	--

7.3 Configure MAIL and FTP target (router systems)

Output	Input	Note
Fig. 33: Mail (name and IP of your SRS) 	Enter hostname and IP address of your Access Server see below table: Select radio-button LAN Click Save Click OK in the pop-up window	

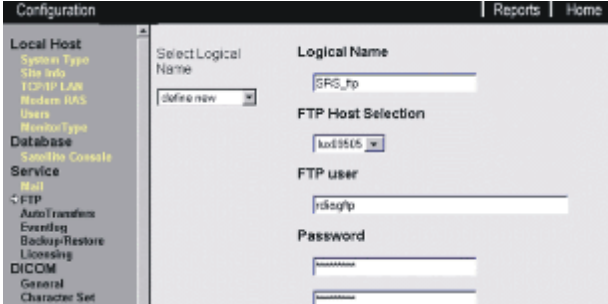
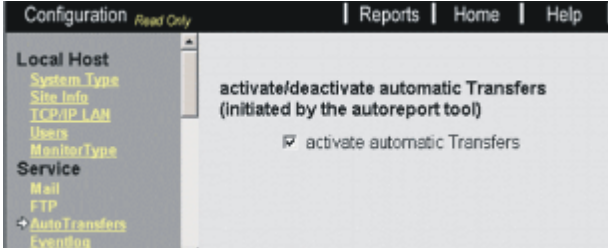
Country	Rdiag server name	IP of Access Server
Europe, Africa	LUX09505	194.138.39.18
North- / South America	SRSACC1	129.73.116.92
Asia	sgpt806x	194.138.243.178

For MR 2004A Rev. B it will be possible to import proactive templates from the Siemens Remote Server for syngo MR systems. You must use the logical names **SRS_mail** and **SRS_ftp** for the mail and ftp targets in order to be able use the preconfigured templates

Output	Input / Note
Fig. 34: Logical name for mail target SRS_mail 	Define mail target Logical Name SRS_mail Select your Access Server in the pull-down menu "Mail Server Host Selection" e-mail address, refer to the table below Click Save Click OK in the pop-up window

If you don't have a e-mail address for your USC, use one of the following address:

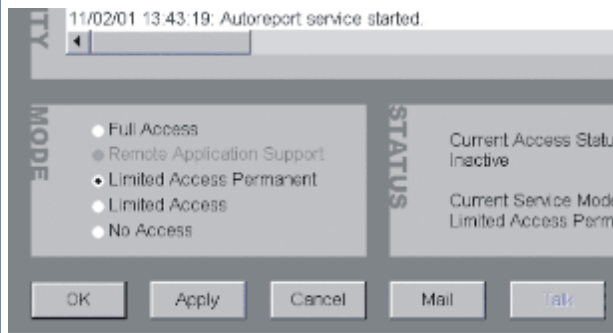
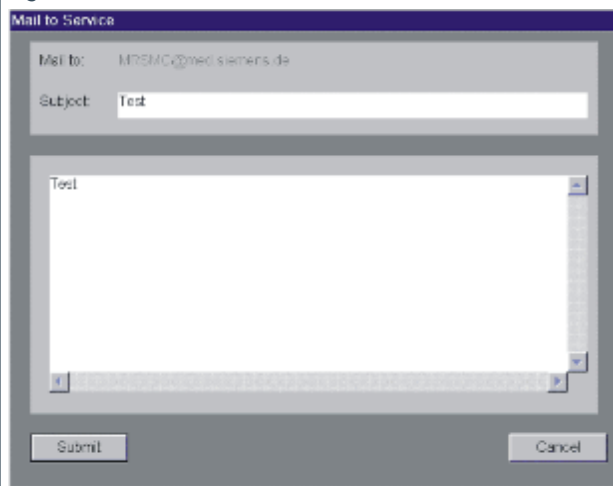
Access Server	e-mail address
Europe, Africa, and Asia	MRSMC.med@siemens.com
North- / South America	MRUPTIME@sms.siemens.com

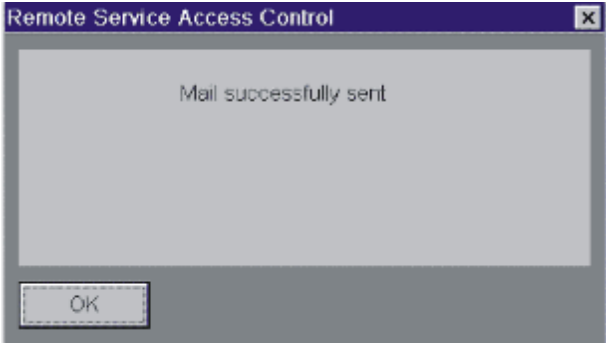
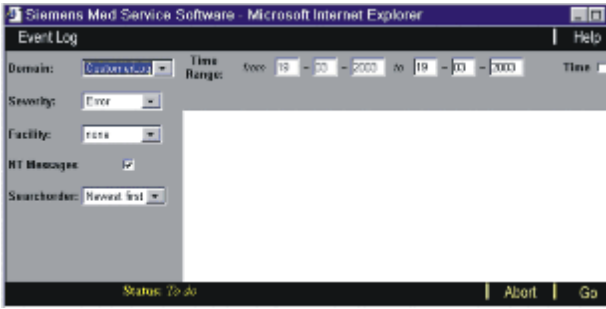
Output	Input	Note
	Skip the first window with >	The host target of mail is used for ftp as well.
<i>Fig. 35: Logical name for ftp target SRS_ftp</i> 	Define ftp target Logical Name SRS_ftp Select your Access Server in the pull-down menu "Ftp Host Selection" FTP user rdiagftp enter the password confirm the password Click Save Click OK in the pop-up window	
<i>Fig. 36: Activate automatic Transfers</i> 	If the customer agrees on remote monitoring Set the check-box: activate automatic Transfers	
For questions regarding the password for the rdiagftp account, please contact your remote diagnostic responsible / USC.		

7.4 Testing connectivity from system to Rdiag server

7.4.1 Send a test mail

- Check at the bottom right that no remote service is active

Output	Input / Note
Fig. 37: Remote Service Access Control 	Select Options, Remote Service Click Mail
Fig. 38: Mail to Service 	Enter subject e.g. SRS test MRxxxxx Text: Please confirm reception of this mail to your.name@siemens.com or under phone number xxx Click Submit

Output	Input / Note
Fig. 39: Mail successfully sent 	Click OK
Fig. 40: Mail successfully sent 	Select Options , Event Log Select Domain: Customer Log

If following error message is displayed

E 19-Mar-03 11:18:35 16 TransferMgr n/a MAIL error.

Receiver: "MRxxx@xxx.siemens.xx"

Filename: "C:\MedCom\service\extract\MRxxxxxx20030319111756Mail.txt".

the customer router is not properly configured for outgoing traffic, please contact the person responsible for the router configuration and the local network administrator.

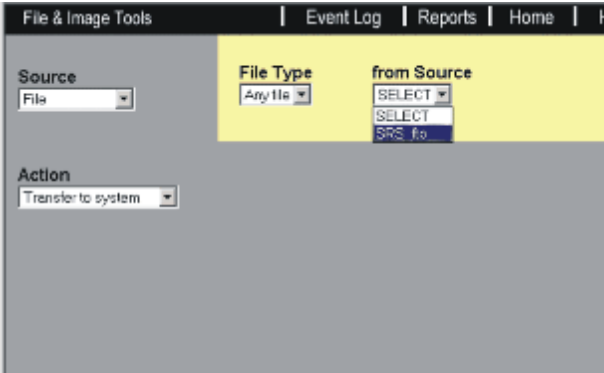
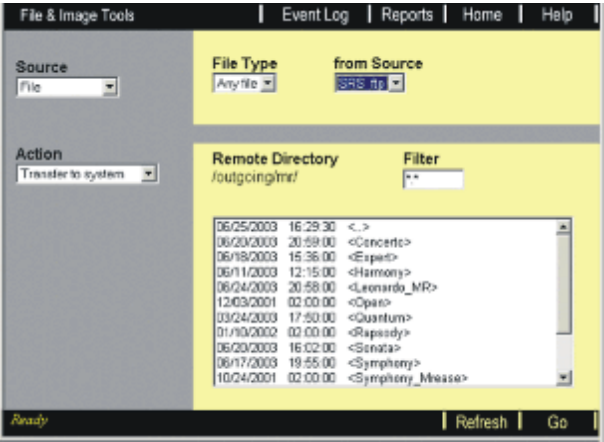
Check to ensure that the configured SIEMENS Remote Server is added to the IP security list, click **Reports**; select **IP Security** below "**SiteInfo**". If the System Management server is missing, refer to (→ Adding host to IP security list / Page 139).

7.4.2 Test FTP login to the Remote Server

7.4.2.1 Prerequisite

The test (→ Send a test mail / Page 104) was successful

- Start the Local Service platform, select File & Image Tools

Output	Input / Note
Fig. 41: File & Image Tools 	Source: select File Action: select Transfer to system File Type: select Any file from Source: select SRS_ftp
Fig. 42: Listing "Remote Directory" 	If a listing of folders is displayed in the frame "Remote Directory" the access to the FTP account of the Remote Server was successful. Click Home Close the Service Platform

- If the test fails, re-enter the password for account "**rdiagftp**" perform the test once more if the test fails again, please contact your Remote Server responsible / USC.

7.5 AutoReport configuration

It is important to transfer data from the site to your USC/TSC/HSC on a regular basis. Quality Assurance and Tune-Up measurements logs, which will be collected in a central database. This allows various kinds of statistic evaluation. The AutoReport function takes over the tedious job of regularly checking for new data and transmitting them manually.



Be aware, that AutoReport can only work correctly, if the following terms are fulfilled:

- LAN available ROUTER / Modem is connected and powered on
- Your site has been entered into the RDIAG server
- Target definition is correct. (Dial mode, escape number, phone number and password (obtained from your USC/TSC/HSC)
- AutoReport is defined and enabled.

The router is required free-of-charge to our customers during the warranty period. It also extends beyond that period if a service contract is concluded.

The planning of the "Remote Diagnostics" environment is performed by the project manager. For detailed information, please refer to the (→ "Planning Guide", <http://cs.med.siemens.de/cms/productinfo/default.asp?startfolder=801281ea3099989879a6fa4ba28d2b0c>).

7.5.1 Configuring QDB report

- **Location of the QDB file:** C:\MedCom\service\html\seso\sitedb\export\export.xml
- **Sample When:** Choose a time, when the system is usually in operation!

The system saves the "AutoReport" configuration, then it sends a test-mail within 5 minutes after completed AR configuration to the mail target destination.

Set-up the Autoreport frequency to daily, to report magnet monitoring data.

- Start the service platform

Output	Input	Note
Fig. 43: 	Click at Auto Report Select TEMPLATE in pull-down menu REPORT	

Output	Input	Note
Fig. 44: Configure qdb auto report (1) Target (2) Name (3) End time (4) Source (5) Filename (6) Mode (7) Sample when (8) Sample when: time (9) Mail (10) FTP target (11) Send when	1. Target: Select "SRS_mail". 2. Name: e.g. qdb 3. End time e.g. 2011 4. Sample source: Select File 5. Filename: Enter location of the "qdb" file, see above path in "Location of the QDB file". 6. Sample what: "Mode:" Select Time Triggered 7. Sample when: Click radio button Daily 8. Sample when: time, see above "Sample When". 9. Send what?, "Mail:" Select as separate file(s) without mails 10. Send what?, "FTP target:" "SRS_ftp" 11. Click radio button At once	
When all of it is configured correctly, click the Save button.		
The system will send a test-mail automatically after saving the AutoReport For results refer to (→ Send a test mail / Page 104).		

7.5.2 Configuring EDB report

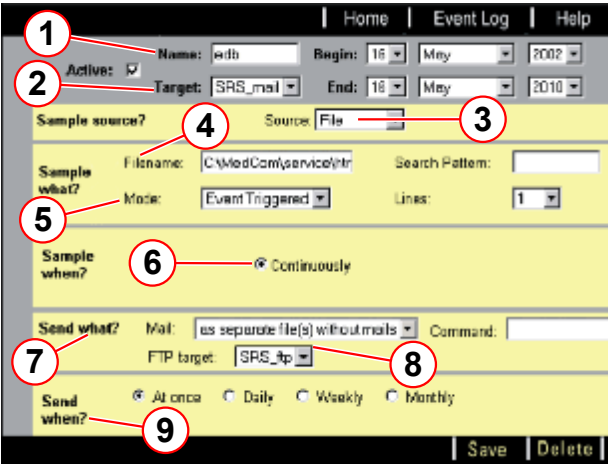
- Location of the EDB file: C:\MedCom\service\html\seso\sitedb\export*.edb

The system saves the "AutoReport" configuration, then it sends a test-mail within 5 minutes after completed AR configuration to the mail target destination.

Set-up the Autoreport frequency to once a week.

- Start the service platform

Output	Input	Note
Fig. 45:	Click at Auto Report Select TEMPLATE in pull-down menu REPORT	

Output	Input	Note
Fig. 46: Configure edb auto report  (1) Name (2) Target (3) Source (4) Filename (5) Mode (6) Sample when (7) Send what (8) FTP target (9) Send when	1. Target: Select "SRS_mail". 2. Name: e.g. edb. 3. Sample source: Select File 4. Filename: Enter location of the "mrs" files, see above path in "Location of the mrs file". 5. Sample what: "Mode:" Select Event Triggered 6. Sample when: Click radio button Continuously 7. Send what?, "Mail:" Select as separate file(s) without mails 8. Send what?, "FTP target:" Select SRSftp 9. Click radio button At once	
When all of it is configured correctly, click the Save button.		
The system will send a test-mail automatically after saving the AutoReport. For results refer to (→ Send a test mail / Page 104).		

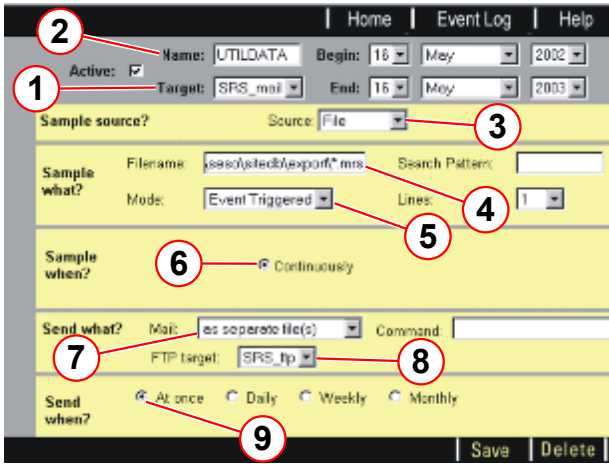
7.5.3 Configure System Utilization reporting

- **Target:** Choose your responsible USC/TSC/HSC for support regarding rdiag server.
- **Location of sysutil files:** C:\MedCom\service\html\sesolsitedb\export*.zip

The system saves the "AutoReport" configuration, then it sends a test-mail within 5 minutes after completed AR configuration to the mail target destination.

- Start the service platform

Output	Input	Note
Fig. 47: 	Click at Auto Report Select TEMPLATE in pull-down menu REPORT	

Output	Input	Note
Fig. 48: Configure System Utilization reporting  (1) Target (2) Name (3) Source (4) Filename (5) Mode (6) Sample when (7) Mail (8) FTP target (9) Send when	1. Target: Select rdiag host, see above "Target". 2. Name: UTILDATA 3. Sample source: Select File 4. Filename: Enter location of the "sysutil files" files, see above path in "Location of the file". 5. Sample what: "Mode:" Select Event Triggered 6. Sample when: Click radio button Continuously 7. Send what?, "Mail:" Select as separate file(s) without mails 8. Send what?, "FTP target:" Select SRSftp 9. Click radio button At once	
When all is correctly configured, click the Save button.		
The system will send a test-mail automatically after saving the AutoReport. For results refer to (→ Send a test mail / Page 104)		

7.6 System Utilization



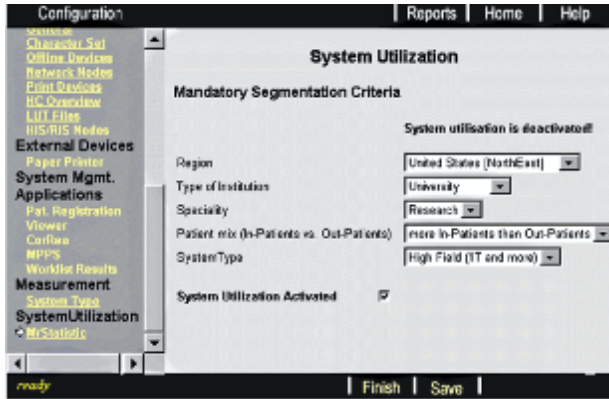
To meet legal aspects, the customer has to sign a contract in advance. Please involve your service organization before you activate System Utilization.

Autoreport for SYSUTIL must be configured before activating System Utilization.



Only activate System Utilization for systems where the customer agrees on "Utilization Management"

The system utilization provides log and report data regarding the usage of syngo MR applications

Output	Input / Note
	Prerequisite to activate System Utilization is that activate automatic Transfers is checked in Configuration, Auto-Transfers.
	Start Configuration in the Service UI, navigate to System Utilization
<i>Fig. 49: Sytem Utilization Activated</i> 	Use the pull-down menu to select region and site profile. If the customer has signed the contract on Utilization Management click the checkbox: System Utilization Activated Click Save to save changes

7.7 System Management Installation

System management is a new remote service for syngo MR systems, it provides following functionality: proactive event reporting and statistical data transfer. In the future it will be used for SW distribution e.g hotfix / buxfix and anti virus patterns.

The installation should be performed at the **MRC** and **MRSC** if available.

7.7.1 Prerequisites

Download the CA test tools from the **CS Intranet**: <http://cs.med.siemens.de>, **For Service, MR, Tools** download the **CA_Com_Check_Tool.exe** to your local service PC.

- Insert a blank floppy into the floppy drive
- Double-click the downloaded file to extract the file by clicking "Unzip". Check if 9 files are successfully extracted (popup message). Click Close to continue.
- Remove the floppy and label it "CA Com Test"

Enable IP security for the CA Unicenter servers in your region

Region: Europe and Africa	server name	IP address
Event Management Server	FTHW9MTA	194.138.39.20
Asset Management Server	FTHW9MUA	194.138.39.21
Software Distribution Server	FTHW9MVA	194.138.39.22

Region: North- / South America and Asia	server name	IP address
Event Management Server	SRSSQL01	129.73.116.90
Asset Management Server	SRSSQL02	129.73.116.93
Software Distribution Server	SRSSQL03	129.73.116.94

Enable the system management servers for your region to the IP security list, use following procedure:

- At the "MRC", simultaneously press **Ctrl - Esc**. Select **Start, Advanced User**; enter the advanced user password; click **OK**
- Simultaneously press **Ctrl - Esc**, click **Start, Run cmd** enter the command: **runas /user:administrator cmd** (enter the administrator password) In the newly started cmd tool you have administrator rights.

for Region: Europe and Africa type following commands in the cmd tool:

```
ipseccmd -f (0+194.138.39.20) ↵
ipseccmd -f (0+194.138.39.21) ↵
ipseccmd -f (0+194.138.39.22) ↵
```


for Region: North- / South America and Asia type following commands in the cmd tool:

```
ipseccmd -f (0+129.73.116.90) ↵
ipseccmd -f (0+129.73.116.93) ↵
ipseccmd -f (0+129.73.116.94) ↵
```

7.7.2 Routing

As the optimal solution, the customer network administrator has configured the route from the default gateway to the SIEMENS server for your continent.

- Skip this step if the SRS router is configured as default gateway. Do NOT add permanent routes if there is no entry under "Persistent Routes" of the logfile `routing.log` under `C:\Medcom\MrisiteData`)

Logfile example **without** any permanent routes:

```
.....
Persistent Routes:
    None
```

Logfile example **with** permanent routes:

```
.....
Persistent Routes:
    Network Address      Netmask                Gateway Address      Metric
    194.138.39.0         255.255.255.255       xxx.xxx.xxx.xxx      1
or
    129.73.116.0         255.255.255.255       xxx.xxx.xxx.xxx      1
or
    194.138.243.0        255.255.255.255       xxx.xxx.xxx.xxx      1
```

The below **example** means **use the customer router / gateway with IP: xx.xx.xx.xx** for the SRS network as destination. E.g. (LUX09505 **194.138.39.18**)
194.138.39.0 route to the whole SRS subnet, and **-p** makes the entry persistent)

Add permanent routes for system in region **Europe and Africa** to the **CA Server**:

```
route add -p 194.138.39.0 mask 255.255.255.0 xx.xx.xx.xx ↵
```

Add permanent routes for system in region **North- / South America and Asia** to the **CA Server**:

```
route add -p 129.73.116.0 mask 255.255.255.0 xx.xx.xx.xx ↵
```

(The Asia server "sgpt806x" automatically routes the IP addresses: 129.7.116.90, 129.7.116.93, and 129.7.116.94) via a VPN tunnel to the CA server)

7.7.3 Check connectivity to the System Management server

- Simultaneously press **Ctrl - Esc**, click **Start, Run cmd**

or use the following script from the floppy "CA Com Test" by entering:

for **Europe, Africa**

a:\CheckCA_Europe.BAT

for **North America, South America, Asia**

a:\CheckCA_US.BAT



Important: Only continue the installation /configuration if the connectivity to the System Management server is successful.

7.7.4 Required ports for System Management (CA Unicenter) Communication

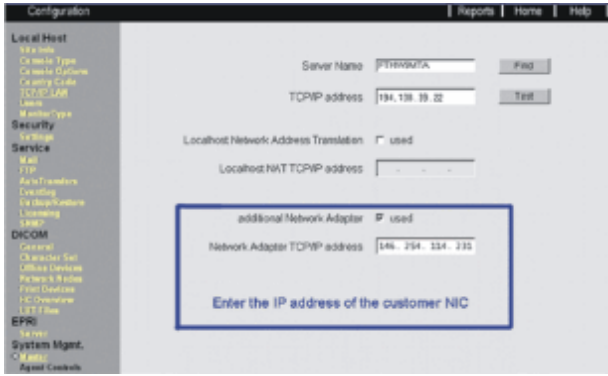
Provide following information to the customer network administrator **if the connectivity test failed** to the System Management server or contact your SRS responsible.

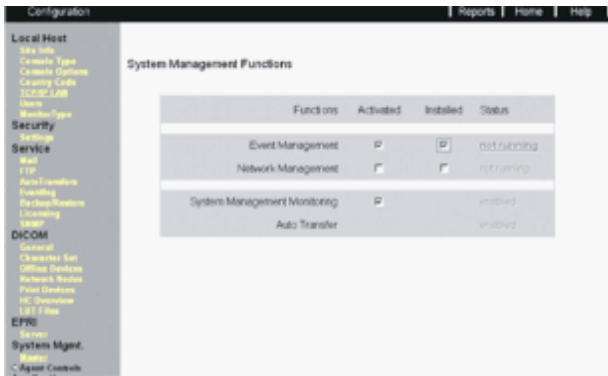
Protocol	Port	Description
TCP	4105	CA Messaging (TCP)
UDP	4104	CA Messaging (UDP)
TCP	8222-8235	Data Transport Agent (DTA)
TCP	7001	Common Communication Interface (CCI)
TCP	1721	CCI Remote Daemon
TCP	4718-4724	Software Delivery (TCP)
UDP	4725	SD Job Check (UDP)

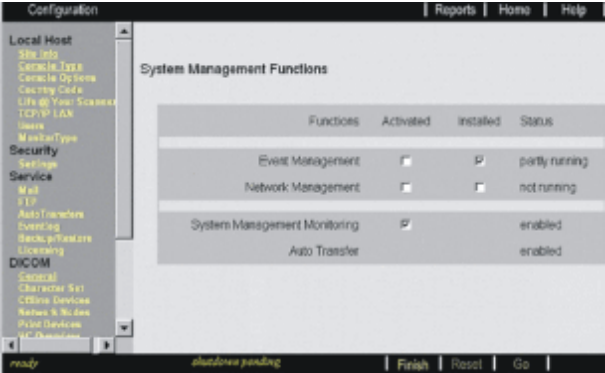
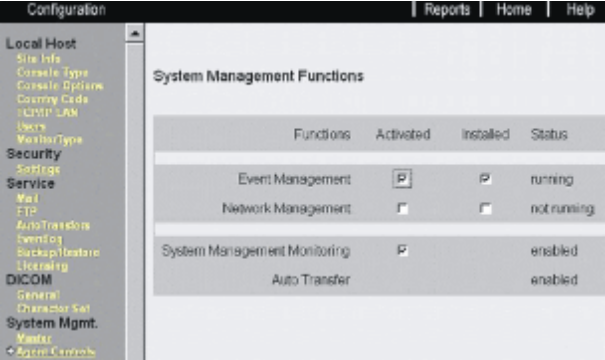
7.7.5 Master

Region: Europe and Africa	server name	IP address
Software Distribution Server	FTHW9MVA	194.138.39.22

Region: North- / South America and Asia	server name	IP address
Software Distribution Server	SRSSQL03	129.73.116.94

Output	Input / Note
<i>Fig. 50: CA Master</i> 	In this mask, you configure the Server Name / IP addresses of the System Management Master for your Region see above table. Click the checkbox "additional Network Adapter" if console type is MRC, enter the TCP/IP address of the network adapter used for SRS (Siemens Remote Services) in the field "Network Adapter TCP/IP address". This is usually the hospital network address of the host.
If NAT (Network Address Translation) is used, click checkbox "Local host Network Address Translation". As "Local NAT TCP/IP address", enter the NAT TCP/IP address assigned in the SRS server system configuration. Ask the local USC for this information.	

Output	Input / Note
<i>Fig. 51: CA Functions</i> 	In this mask, you may stop and start the System Management agents. Functions Select Event Management and System Management Monitoring). Click the 'GO' button to start the Sytem Management Installation/configuration.
The "System Management" installation / configuration takes ~25 minutes, please wait.	

Output	Input / Note
Fig. 52: CA Check if installation was successful 	The installation was successful, when the status has changed to partly running and the message "shutdown pending" is displayed in the in status bar of the Service UI.
System management installation / configuration pop-up will be displayed at the end of the installation.	
Fig. 53: CA Check if installation was successful 	After the system re-boot: Check that: ■ Event Management has the status running ■ System Management Monitoring has the status enabled

System Management Functions

- i

If only one or two checkboxes are marked, CA Unicenter (System Management) is only partly installed. In this case, deinstallation of CA Unicenter (System Management) is required!
- i

If Auto Transfer is disabled and Event Management is selected, the collected data will not be transferred to the System Management unless a manual transfer is performed. This may result in data loss since the buffer size is limited (ring buffer).

7.8 Troubleshooting and useful information

7.8.1 System Management Installation failed

In case that the System Management Installation fails perform a manual de-installation as described below:

- Login as user administrator
- click **Start, Run** cmd
- **cd C:\MedCom\service\mwtools**
- **sysmgmt_install.exe -r**

- Check the status of CA Unicenter (System Management) by using the following command in the cmd tool window:

```
sysmgmt_install.exe -c<ENTER>
```

The following shows a successful deinstallation of CA UniCA Unicenter (System Management)

```
C:\MedCom\service\mwtools>sysmgmt_install.exe -c
```

```
*****
```

```
* Systemmanagement Managed Node Package - Version VE10D
```

```
*****
```

```
OS=WinXP/.Net
```

```
-----
```

```
Check Installation-Account
```

```
Account=administrator
```

```
-----
```

```
Check Registry-Keys
```

```
All Registry-Keys are present
```

```
-----
```

```
addr #0 is [192.168.1.1]
```

```
addr #1 is [192.168.2.1]
```

```
*****
```

```
* Systemmanagement Managed Node Package - Version VE10D Status:
```

```
*****
```

```
Backup Management: not installed
```

```
Event Management: not installed
```

```
Asset Agent not installed
```

```
Event Agent not installed
```

```
Hardware Agent not installed
```

```
Log Agent not installed
```

```
Software Delivery Agent not installed
```

```
Software Transport Agent not installed
```

```
System Agent not installed
```

```
Performance Management: not installed
```

```
Network Management: not installed
```

```
Rules/Filters Version: not installed
```

```
C:\MedCom\service\mwtools>
```

- If any entry from the "System Management Managed Node Package" shows installed, the program "sysmgmt_install.exe -r" must be run again. The status of CA Unicenter (System Management) has to be checked as well after the reboot.

If it still fails, contact the SRS help desk (Tel. xx49-9191-188080-191)

For problems within the US, call the US Uptime Service Center / technical support group.

For further investigation get the following logfile C:\sysmgmt\EvtMgmtI-NIT.log and C:\TEMP\sysmgmt_install.log.

8.1 Re installation of software fails

If the SW re installation fails e.g. due to a interrupt during the initial installation part of CD 1/2. The install software was partly portated to the system disk. At the next re installation, the install SW supposes that the SW installation is already is performed with a dead lock as result.

8.1.1 Formatting the system disk

Work around:

- Insert the CD entitled "syngo MR CD 1/2" and reboot the system. Press "Alt-Ctrl-Del" simultaneously
- When the screen with the error message is displayed, press "CTRL-C" to cancel the installation.

Output	Input	Note
The DOS-Prompt appears on the left side of the screen.	a: ↵	If the drive letter is "C:\", then change to drive A:\ MS DOS including FDISK is loaded on RAM during the boot from CD 1/2. The ram content simulates a DOS floppy . No floppy is needed!
A:\	cd utils ↵	
List of options ...	fdisk ↵	Starts the partitioning program.

- Select the 4-th menu item, to get an overview of the existing partitions. Check now the file system type of the first partition.
- Press "**ESC**" to leave the current menu.
- Choose menu item "**3**" to delete Partitions.
- Depending on partition type, choose one of the items 1 to 4 to delete your first partition. In case you typed in a false menu item, "ESC" will bring you back to Step 8. This also applies, if message "NO PARTITION TO DELETE" appears after choosing a menu item.
- When you choose the right type for your first partition, the warning "DATA IN THE ... PARTITION WILL BE LOST" and "WHAT ...PARTITION DO YOU WANT TO DELETE [1]" appears
- Press "ENTER (↵)".

If a request for the partition name is asked, type in "MED_SYSTEM", which will also be displayed at the top of the screen.

- Continue with "Y" for "Yes"



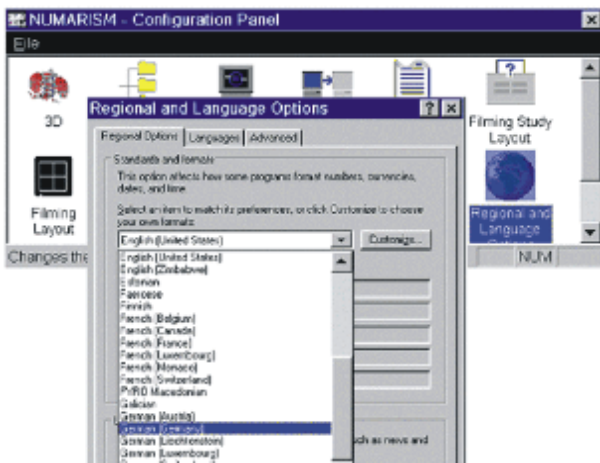
For some keyboard layouts the "Z" key will produce the "Y" in the DOS mode.

The message " ...PARTITION DELETED" appears.

- Exit the fdisk program by pressing "ESC" twice, press any key to reboot.

Try again to re install software, make sure that the "syngo MR CD 1/6" is still inserted in the CD-ROM drive.

8.2 System language and weight

Output	Input	Note
Fig. 54: 	Click Options Click Configuration Click Regional and Language Settings The "Regional Options" tab must be selected Select system language in the pull-down menu. Click Apply, Click OK	The system language must be set via the customer configuration platform within the MR application Patient weight unit for the "Patient Registration" input field, [lb] or [kg] is in relation to selected system language. English US = 1ball others = kg) Don't set the checkbox: "Set as system default locale" The language settings will be activated after following application restart
	Reboot the system	

8.3 Keyboard

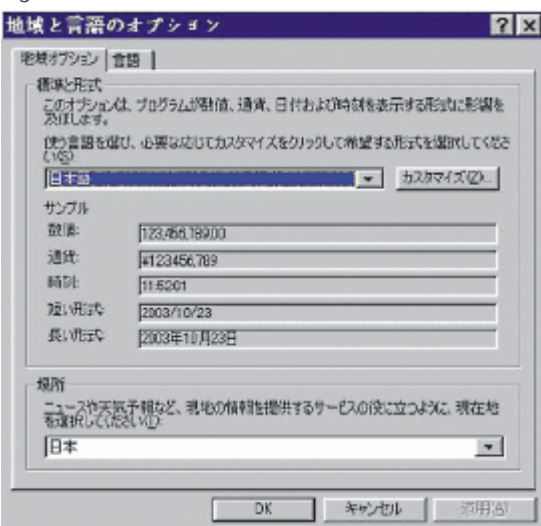
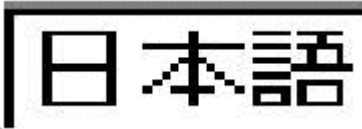
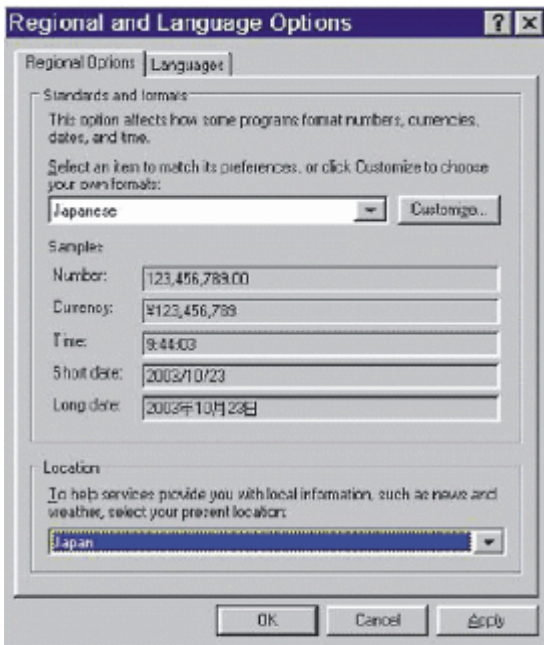
Refer to (→ Keyboard settings and Digit grouping symbol for meduser / Page 63).

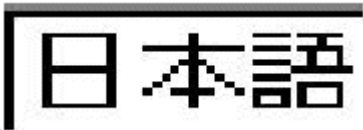
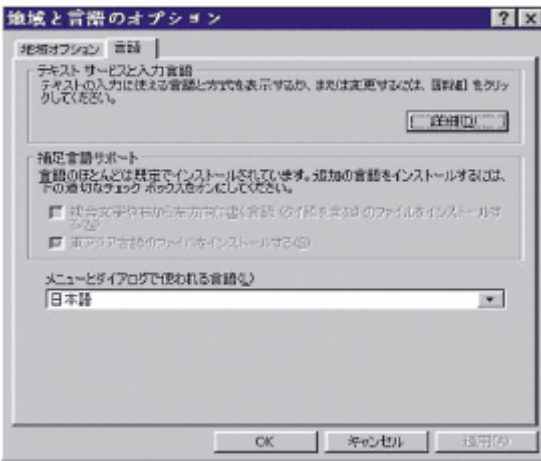
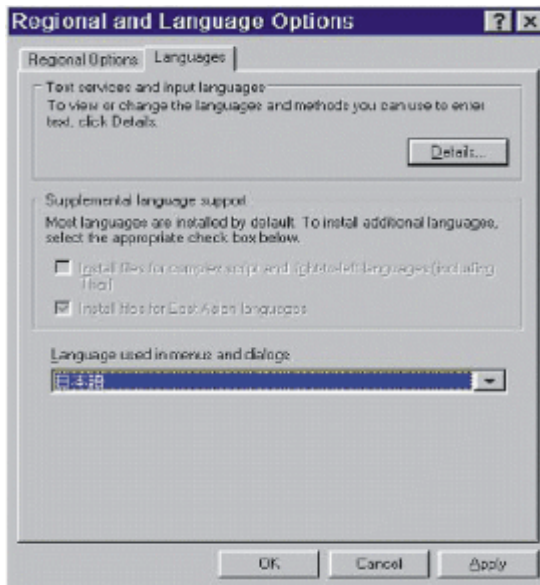
8.4 Language Settings Japanese

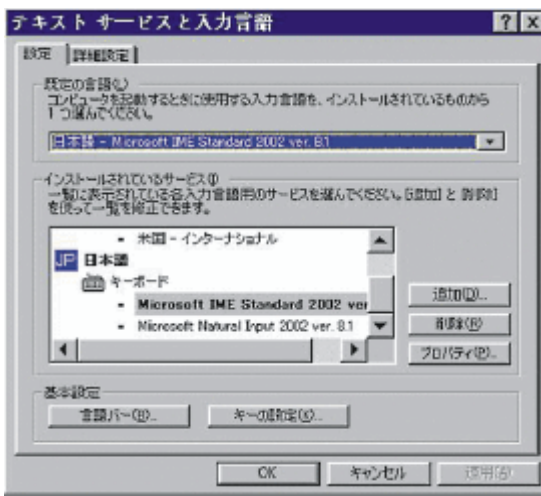
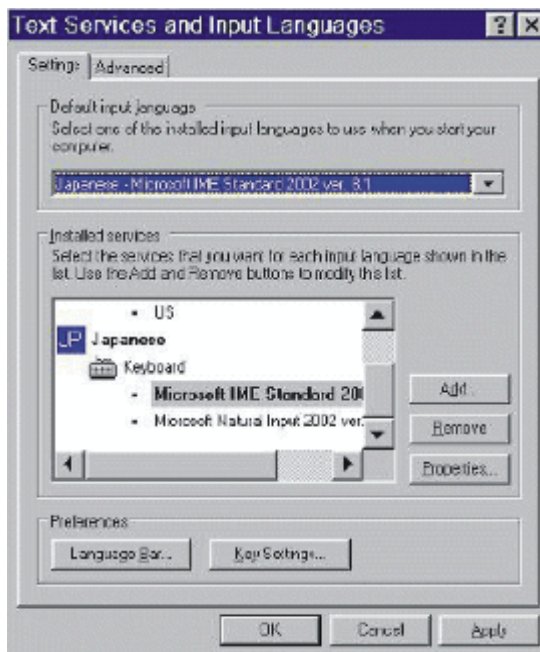
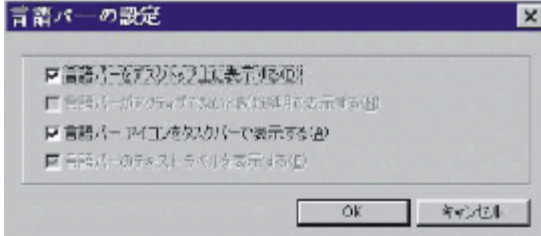
8.4.1 Switch language from international (english-US) to japanese

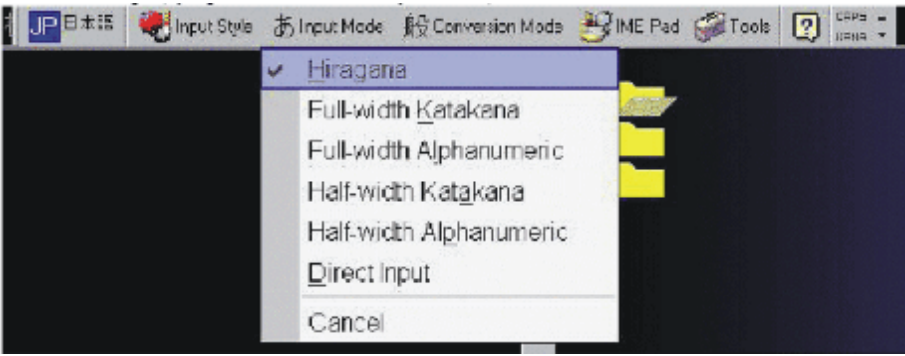


The following steps can only be done, if the japanese version of WinXP (syngoMR 2004A Customer CD6/6) was used for the installation.

ACTION		MRC	MRSC	
■ Select Options -> Configuration in the main menu of the Meduser platform		☐	☐	
■ Double-click the Regional and Language Settings icon		☐	☐	
■ In the highlighted field choose Output popup in Japanese, if it was already switched during Windows XP Installation. Fig. 55:  the first Item (see below)  (Hint: scroll up, if the following symbol is not present)		Else the international view of the popup. Fig. 56: 	☐	☐

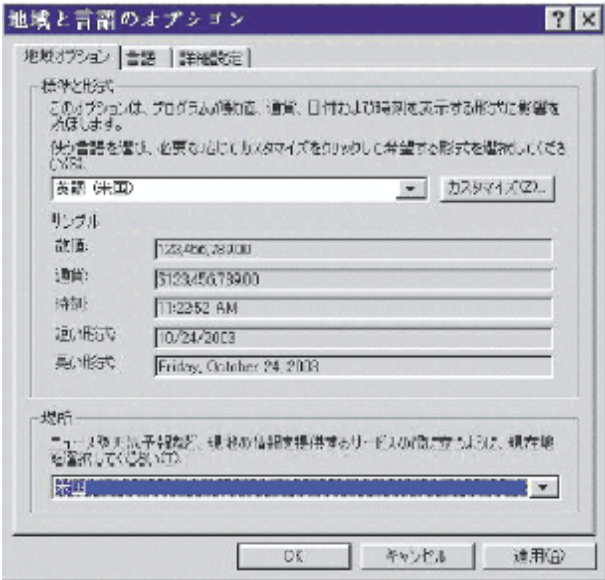
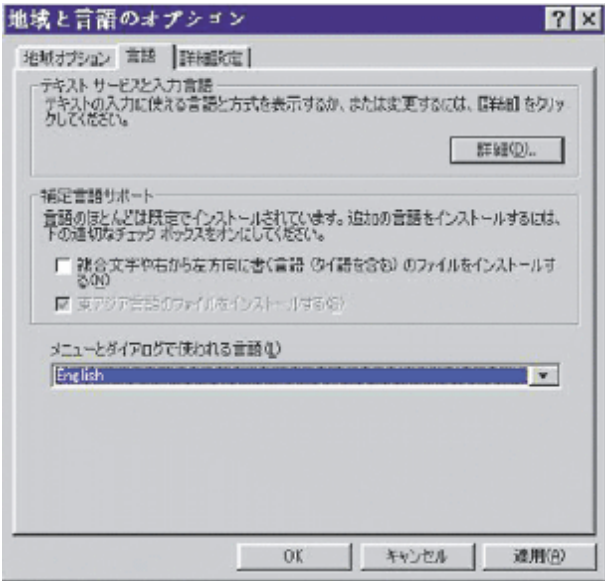
ACTION		MRC	MRSC
Note: Be sure to use the long format for the date. It must look like yyyy/mm/dd or yyyy-mm-dd. If it ist not set you have to change it. For it press the Customize button and select the correct setting for the date. Switch to the 2. tabcard "language" and choose the same Japanese symbol as in the first tabcard. Japanese symbol: 			
Japanese view Fig. 57: 	or international (english-US) view Fig. 58: 	☐	☐
■ Click on the "Details" Button and make following settings for your Keyboard:			

ACTION		MRC	MRSC
Japanese view Fig. 59: 	or international (english-US) view Fig. 60: 	☐	☐
Click the Language Bar Button under preferences and select the two checkboxes like shown below:			
Japanese view Fig. 61:  Press the OK button.	or international (english-US) view Fig. 62:  Press the OK button.	☐	☐
First press OK, next Apply. Note: Pressing "OK" in the "Regional and Language Options" Tab Card starts an automatic reboot, wait two steps with this action.		☐	☐
Open the Windows Explorer (It could be necessary to switch to "Advanced User" to open an Explorer)		☐	☐
Delete the file C:\MedCom\config\Mri\MRSystem\MrVCheck.ref		☐	☐
Press the OK Button and accept all Messages, which will be shown with "YES" or "OK". A Reboot follows!		☐	☐
Login as meduser (if autologin is not enabled) after the reboot and wait until the main UI is up.		☐	☐

ACTION	MRC	MRSC
IME settings: (by right mouse-click on "Input Mode"). Fig. 63: 	☐	☐
- Be sure that in the Configuration (Local service) the Dicom Character Set is set to "ISO 2002 IR13" for Japanese! - Start local service - Go in Configuration - Select under DICOM the link Character Set - Choose now the radio button "allow multiple character sets" - Select the Character set "ISO 2002 IR13" - Press "Save" and confirm the popup. - Restart the system by pressing "Home". Confirm the restart popups.	☐	☐

8.4.2 Switch language from japanese to international (english-US)

ACTION	MRC	MRSC
Select Options -> Configuration in the main menu of the Meduser platform	☐	☐
Double-click the Regional and Language Settings icon	☐	☐

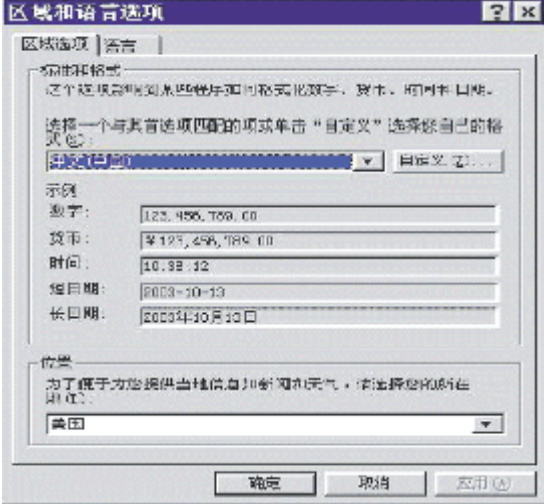
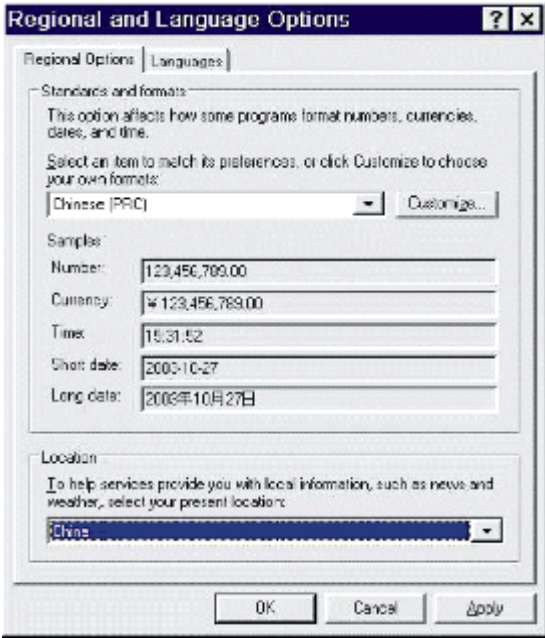
ACTION	MRC	MRSC
Change on the first tabcard (Regional Options) the settings like in the picture below English - US. Fig. 64: 	☐	☐
Change on the second tabcard (Language) the settings like in the picture below English – US. Fig. 65: 	☐	☐
First press OK, next Apply. Hint: Pressing "OK" in the "Regional and Language Options" Tab Card starts an automatic reboot, wait two steps with this action.	☐	☐
Open the Windows Explorer (It could be necessary to switch to "Advanced User" to open an Explorer)	☐	☐

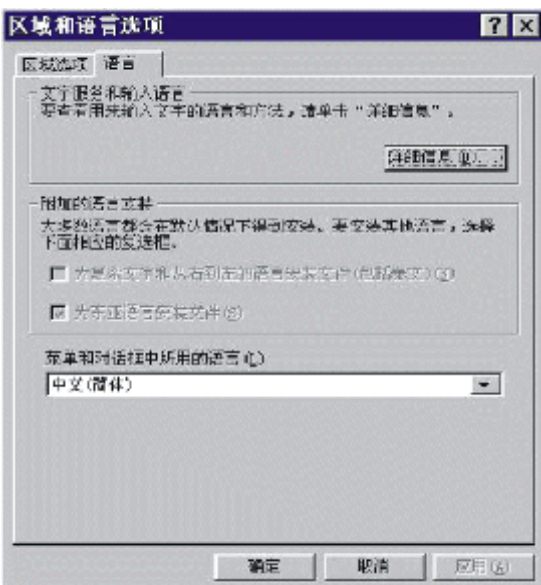
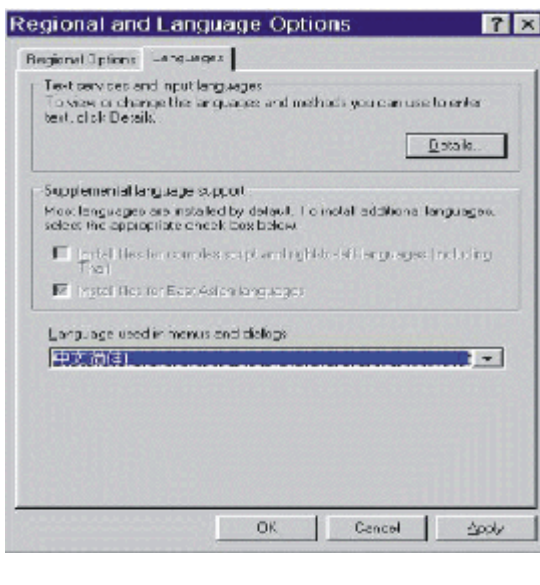
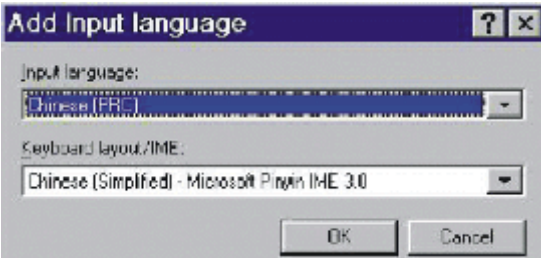
ACTION	MRC	MRSC
■ Delete the file C:\MedCom\config\Mri\MRSystem\MrVCheck.ref	☐	☐
■ Press the OK Button and accept all Messages, which will be shown with "YES" or "OK". A Reboot follows!	☐	☐

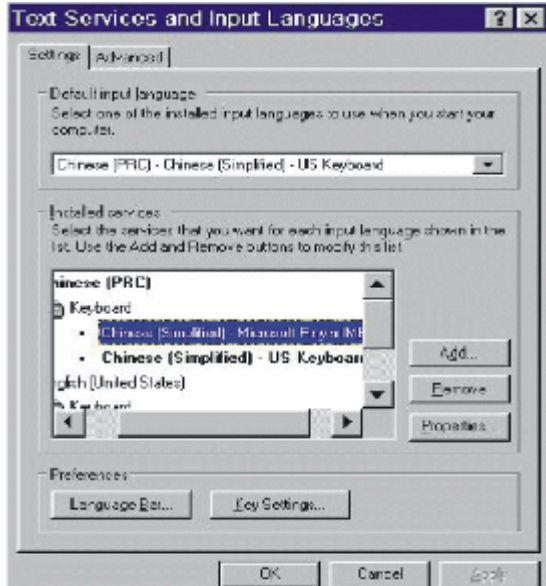
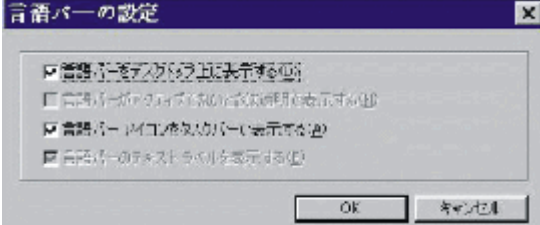
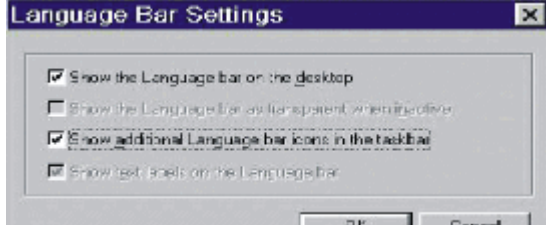
8.5 Language Settings Chinese

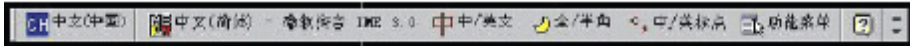
8.5.1 Switch language from international (english-US) to chinese

Note: The following steps can only be done, if the chinese version of WinXP (syngoMR 2004A Customer CD5/6) was used for the installation

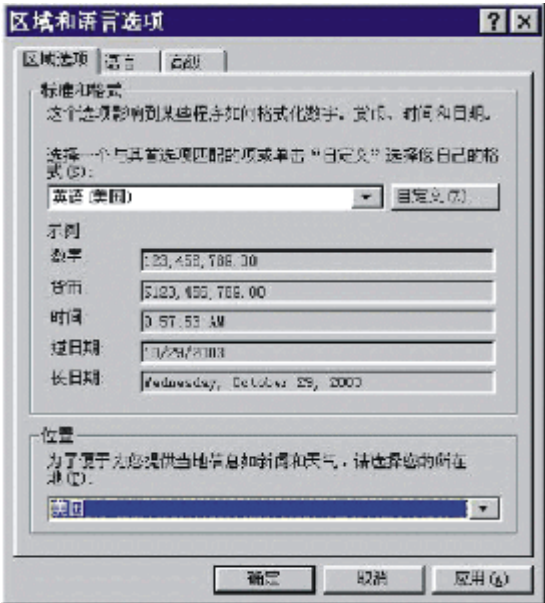
ACTION		MRC	MRSC
■ Select Options -> Configuration in the main menu of the Meduser platform		☐	☐
■ Double-click the Regional and Language Settings icon		☐	☐
■ In the highlighted field choose Output popup in Chinese, if it was already switched during Windows XP Installation. Fig. 66: 	Else the international (english-US) view of the popup. Select the below shown entries. Fig. 67: 	☐	☐
the following Item (see below)			
中文(中国) (Hint: scroll up, if the following symbol is not present)			
■ Switch to the 2. tabcard "language" and choose the same Chinese symbol as in the first tabcard. Chinese symbol: 中文(中国)		☐	☐

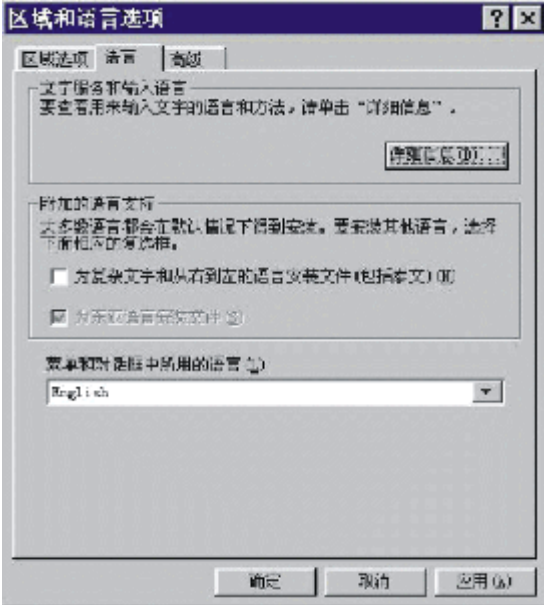
ACTION		MRC	MRSC
Chinese view Fig. 68: 	or international (english-US) view Fig. 69: 		
- Click on the "Details" Button Note: If the necessary Keyboard Layout/IME is not present, press the "Add" Button and select (see following popup) and press OK. Fig. 70:  - Make than the following settings for your Keyboard:		☐	☐

ACTION		MRC	MRSC
Chinese view Fig. 71: 	or international (english-US) view Fig. 72: 		
Click the Language Bar Button under preferences and select the two checkboxes like shown below:		☐	☐
Japanese view Fig. 73:  Press the OK button.	or international view Fig. 74:  Press the OK button.		
First press OK, next Apply. Note: Pressing "OK" in the "Regional and Language Options" Tab Card starts an automatic reboot, wait two steps with this action.		☐	☐
Open the Windows Explorer (It could be necessary to switch to "Advanced User" to open an Explorer)		☐	☐
Delete the file C:\MedCom\config\Mri\MRSystem\MrVCheck.ref		☐	☐
Press the OK Button and accept all Messages, which will be shown with "YES" or "OK". A Reboot follows!		☐	☐

ACTION	MRC	MRSC
Login as meduser (if autologin is not enabled) after the reboot and wait until the main UI is up.	☐	☐
You can switch the output language of IME to Chinese by a mouse-click on it: (See the red symbol in the middle). Fig. 75:  or to English Fig. 76: 	☐	☐

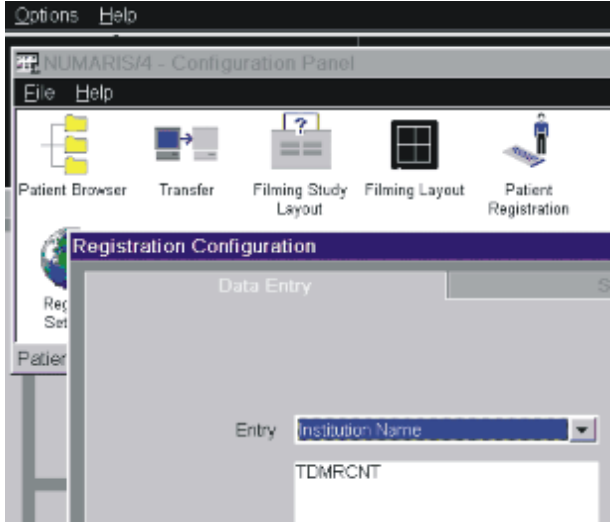
8.5.2 Switch language from chinese to international (english-US)

ACTION	MRC	MRSC
Select Options -> Configuration in the main menu of the Meduser platform	☐	☐
Double-click the Regional and Language Settings icon	☐	☐
Change on the first tabcard (Regional Options) the settings like in the picture below English - US. Fig. 77: 	☐	☐

ACTION	MRC	MRSC
■ Change on the second tabcard (Language) the settings like in the picture below English – US. Fig. 78:	☐	☐
■ First press OK, next Apply. Hint: Pressing “OK” in the “Regional and Language Options” Tab Card starts an automatic reboot, wait two steps with this action.	☐	☐
■ Open the Windows Explorer (It could be necessary to switch to “Advanced User” to open an Explorer)	☐	☐
■ Delete the file C:\MedCom\config\Mri\MRSystem\MrVCheck.ref	☐	☐
■ Press the OK Button and accept all Messages, which will be shown with “YES” or “OK”. A Reboot follows!	☐	☐

8.6 Institution name

The Customer UI must be up and running.

Output	Input	Note
Fig. 79: 	Click Options Click Configuration Click Patient Registration Click Institution Name	The Institution name will be part of the DICOM header

8.7 Regional settings

8.7.1 Date and Time syntax

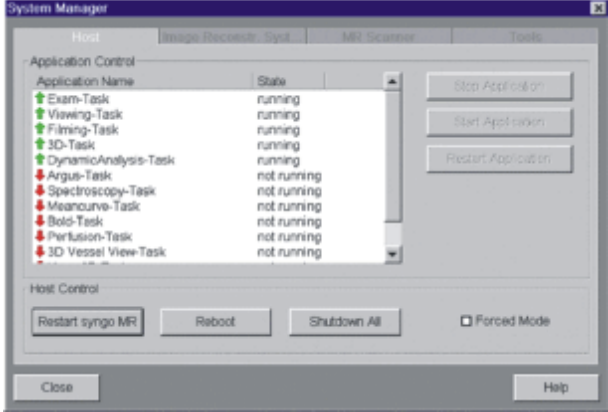
To configure the date and time syntax, proceed as follows: select Options, Configuration, Regional Settings, **Time or Date** in the MR application. It is recommended to select a 4 digit syntax for the year in the "Short date style" pull-down menu to ensure the correct patient age on film printouts (e.g., for patients born in the 19th century or infants).

When you are finished defining the regional settings, restart the MR system to activate the changes.

8.8 System Manager

8.8.1 Restart Application

In case the MR application has crashed or is frozen, try starting the System Manager as follows:

Output	Input	Note
	Simultaneously press [Ctrl]-[Alt]-s	
Fig. 80: 	Click Restart syngo MR	

8.9 Copy SIEMENS Programs to the customer tree

Use the following procedure to copy SIEMENS applications to the customer tree. The customer can modify these templates according to his own requirements.

- Start the Service platform

Output	Input	Note
	cmd.exe	
	InstallProtocols	



If there are existing protocols in the customer tree, the "InstallProtocols" command creates a new tree without overwriting the existing one. The new SIEMENS program folders are numbered "001", "002", etc.

8.10 IP Security

The syngo-MR systems are protected with a service "ipsec".

Other network nodes, like DICOM nodes, can only access the DICOM services that your console provides if they are included in the "Generic Transport Filter" of ipsec.

To list nodes that are allowed to access your console, click **Reports**; select **IP Security** below **SiteInfo**.

8.10.1 Adding host to IP security list

To add a node manually into the IP security list, use following procedure:

At the "MRC", simultaneously press **Ctrl - Esc**. Select **Start, Advanced User**; enter the advanced user password; click **OK**

Simultaneously press **Ctrl - Esc**, click **Start, Run cmd** enter the command: **runas / user:administrator cmd** enter the password

In the newly started cmd tool you have administrator rights.

(replace the address "146.254.222.333" in the following **example** with the IP address of your remote node)

Input in the cmd window	Output
ipseccmd -f (0+146.254.222.333)	The command completed successfully

8.10.2 Deleting all IP security settings

To delete all nodes manually from the IP security list, use following procedure:

At the "MRC", simultaneously press **Ctrl - Esc**. Select **Start, Advanced User**; enter the advanced user password; click **OK**

Simultaneously press **Ctrl - Esc**, click **Start, Run cmd** enter the command: **runas / user:administrator cmd** enter the password

Input in the cmd window	Output
ipseccmd -u	The command completed successfully

8.10.3 Adding all configured hosts to IP security list

To add configured nodes manually into the IP security list, use following procedure:

At the "MRC", simultaneously press **Ctrl - Esc**. Select **Start, Advanced User**; enter the advanced user password; click **OK**

Simultaneously press **Ctrl - Esc**, click **Start, Run cmd** enter the command: **runas / user:administrator cmd** enter the password

Input in the cmd window	Output
C:\MedCom\utils\SetIPSec.bat /install	... Setting dynamic security policy for details see C:\MedCom\log\SetIPSec.log

8.11 Context Vision

The Context Vision filter option requires an additional license key that is defined by the **dongle** and the system **serial number**.

- Start the service platform (service level 7 is required)

8.11.1 Requested data for cvlicense generation

Output	Input	Note
	Click Utilities Select Escape to OS	
<i>Fig. 81: Get Host ID</i> 	Type χωλιχενσε ↵	

Your notes:

Product	
Serial n/o	
Host Name:	
Dongle ID	Without 7- , 9- or Dallas dongles without 8-

All licenses are handled by TD HSC 1 in Germany (e-mail: TDHSC@med.siemens.de, phone: **+49-9191-188080**). (Please ask for the "**call center**".)

8.11.2 Activate license for Context Vision

Output	Input	Note
	Click Utilities Select Escape to OS	
<i>Fig. 82: Enter license for Context Vision</i>  Enter license for Context Vision	Type cvlicense [license key] ↵	

What to do next?		
If you are configuring the MRC or MRSC, refer to (→ Restoring site specific data MR 2004A Rev. B on MR 2004A Rev. B / Page 84) 	If you are installing / adding an MRSC, refer to (→ Import configuration data from CD-R to MRSC / Page 61) 	If you are performing the first configuration, refer to (→ Configuration / Page 89) 

8.12 MR auto login

8.12.1 Access NT login

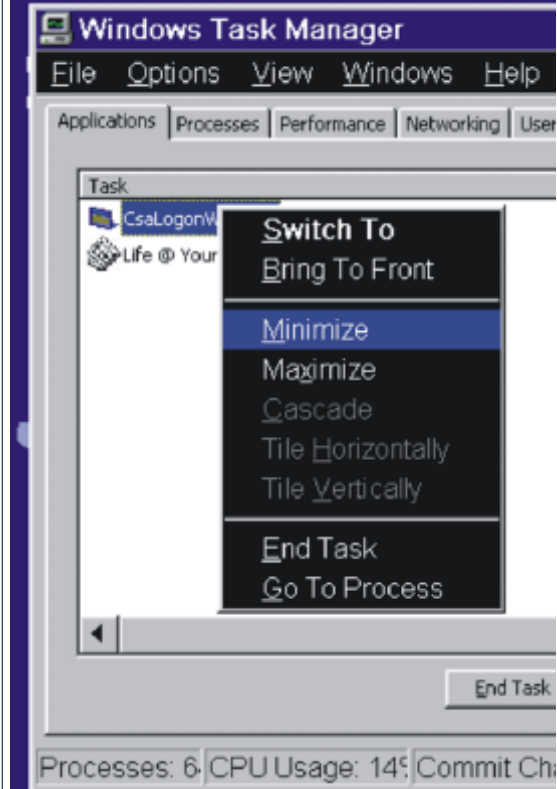
Press SHIFT during the "blue screen" boot phase to access the NT login.



To disable/enable MR autologin, refer to On-line Help texts "Configuration", "Users".

8.13 syngo startup / shutdown counter

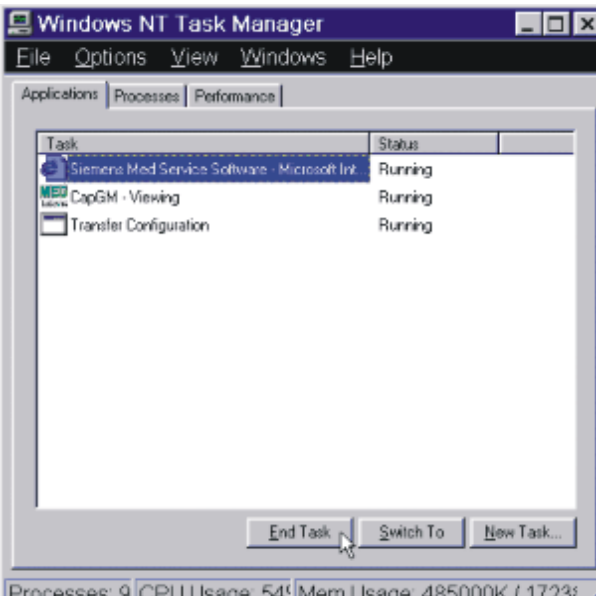
8.13.1 How to terminate the syngo counter

Output	Input	Note
Fig. 83: syngo counter 	Terminating the syngo counter: Press Simultaneously Ctrl, Shift, and ESC Select the tab Applications right click CsaLogonWindow select Minimize	To view actions during startup and shutdown, you can minimize the syngo counter

8.14 Terminate deadlock of crashed Service UI

Internet Explorer can cause a system deadlock during a service session.

8.14.1 How to terminate the Service IU (Internet Explorer)

Output	Input	Note
Fig. 84: Windows NT Task Manager  The screenshot shows the Windows NT Task Manager window with the 'Applications' tab selected. The task list contains three entries: 'Siemens Med Service Software - Microsoft Int.' (Running), 'CapGM - Viewing' (Running), and 'Transfer Configuration' (Running). The 'End Task' button at the bottom is highlighted with a mouse cursor.	Simultaneously press "Ctrl-Alt-Delete" Click the "Applications" tab Select the Siemens Med Service Software task in the task list Click "End Task"	The service platform crashes, the system screen is white (The User Interface will recover within one minute)

8.15 Expert - i

Following is only valid for systems with the "**Expert-i**" option To check if the "MR_Expert_i" licence is available, use the Service: UI, Configuration, Service, Licensing.

- Log in as **administrator** at the **MRC or MRSC**

Output	Input	Note
	Insert the "installtion CD" for the option "Expert-i" to the CD or DVD ROM-drive	
	Open the "Windows file Explorer"	
	Navigate to CD or DVD drive (R:)	
	Double-click the Host.exe	This installs the Expert-i option
	Double-click the Host.exe	Due to the MR firewall the skript must be run once more
	Re-boot, login as meduser	The Expert-i icon will be found in the status-bar

8.16 Measurement system configuration matrix

SW Product No.:	Product Name
024	Harmony
087L	Harmony Impact (Large cabine)
087S	Harmony Impact (Small cabine)
047	Symphony
089	Symphony Vision
085	Sonata
086	Sonata Vision
007	Concerto
083	Allegra
091	Trio

SW Product No.:	024	087L	087S	047	089	085	086	007	083	091
MPCU										
PPC/266 MHz	x	x	x	x	x	x	x	x	-	-
PPC/300 Mhz	x	x	x	x	x	x	x	x	x	x
RX channels										
RX4	x	x	x	x	x	x	x	x	x	-
2. RX4	o	o	o	o	o	o	o		-	x
Receiver										
Receiver 1	x	x	x	x	x	x	x	x	x	-
Receiver 2	o	x	x	o	x	x	x	o	o	-
Receiver 3	o	o	o	o	o	o	o		-	-
Receiver 4	o	o	o	o	o	o	o		-	x
Transmitter										
TRA1	x	x	x	x	x	x	x	x	x	x
TRA2				o	o	o	o	-	-	o
EFI										
No EFI_DAS	-	-	-	-	-	-	-	x	-	-
EFI_DAS	-	-	-	-	-	-		x	-	-

SW Product No.:	024	087L	087S	047	089	085	086	007	083	091
EFI_DAS with Sensor	-	-	-	-	-	-	-	X	-	-
Mains frequency										
50 Hz	-	-	-	-	-	-	-	X	-	-
60 Hz	-	-	-	-	-	-	-	X	-	-
Railway frequency (only valid for sites located close to a railway track)										
16 2/3 Hz	-	-	-	-	-	-	-	X	-	-
25 Hz	-	-	-	-	-	-	-	X	-	-
Gradient amplifier										
Platform/400V	X			X						
Platform/800V	X			X						
Cascade/500V	X			X						
Cascade/1000V 300A	X			X						
Cascade/1000V 380A	X			X						
Cascade/1500V 380A	X			X						
Cascade/2000V 380A	X	X	X	X	X					
Cascade/2000V 500A						X	X		X	X
GPA007								X		
Gradient coils										
AS39S	X			X						
AS39SR	X			X						
AS39T	X	X	X	X	X					
AS Sonata						X				
AS Sonata Vision							X			X
GC018								X		
AC 44	-							X	X	-
Stimulation monitor										
GPA_INT	X	X	X	X	X	X	X	-	X	X
dB/dt Look Ahead Guideline										
SAFE model	X	X	X	X	X	X	X		X	X
FDA	X	X	X	X	X	X	X		X	X

SW Product No.:	024	087L	087S	047	089	085	086	007	083	091
IEC	x	x	x	x	x	x	x		x	x
MHW	x	x	x	x	x	x	x		x	x
RF Amplifier										
DORA (RFPK2204)	x	x	x	x	x	x	x			
CORA (RFPK2137)				o	o	o	o	x		
RFPA (LPPA 13040)				o	o	o	o	-	x	-
RFPA (LPPA 13080)				-	-	-	-	-	-	x
RF Resonator										
BC LP rx/tx	x	x	x						-	-
BC CP rx/tx				x	x	x	x		-	-
BC 2 CP rx/tx				x	x	x	x		-	-
BC LP -/tx								x	-	-
RF safety watchdog										
RFSWD	x	x	x	x	x	x	x		x	x
Number of Body Tune Boxes										
One BTB installed	x	x	x						-	-
Two BTBs installed				x	x	x	x		-	-
BTB Types										
BTB 024	x								-	-
BTB 024S			x						-	-
BTB 024L		x							-	-
BTB 047 AS<5				x					-	-
BTB 047 AS>=5				x					-	-
BTB 047 Syngo				x	x				x	-
BTB Sonata						x			-	-
BTB Sonata Vision							x		-	-
BTB 063								x	-	-
BTB 047-2				x	x	x	x		-	-
BTB 3TWB									-	x
RFCI										

SW Product No.:	024	087L	087S	047	089	085	086	007	083	091
CAPTICON								x		
2.DYSCON	x	x	x	x	x	x	x		x	x
Receiver										
LCCS	x	x	x	x	x	x	x		x	x
COSIMUX								x		
Magnet										
OR41		x	x							
OR42		x	x							
OR35					x		x			
OR24	x									
OR70				x		x				
OR64										x
OR92									x	
Sumitomo								x		
Magnet Supervision										
MSU	x	x	x	x	x	x	x		x	x
TCU								x		
Active Shim										
Additional shim	o	o	o	o	o	x	x		x	x
array shim / standard	x	x		x	x	x	x		-	-
array shim / small cabin			x						-	-
mechanic shim								x	-	-
Front end	o	x	x	o	x	x	x	o	-	-
Patient table										
Distance Standard	x			x		x			-	
Distance Impact		x	x						-	
Distance Vision					x		x		-	x
Open-Table								x		
OP Table				o				o		
OHF Table										

SW Product No.:	024	087L	087S	047	089	085	086	007	083	091
LCD-Display+Banane	x	x	x	x	x	x	x			x
DTP								x		
Spectroscopy (structure related to appl. packages)										
Proton spectroscopy				o	o	o	o		x	x
SAMI-Ansteuerung				o	o	o	o			
2.RFPA				o	o	o	o			
Cooling system										
RCA 1 cooler	x	x	x	x	x		x		-	-
RCA 2 cooler						x			x	x
GRACCO								x		
Mobile										
Mobile enforce	o			o		o			-	-
Coils										
CP Flex large	o	o	o	o	o	o	o			
CP Flex small	o	o	o	o	o	o	o			
Flex loop large	o	o	o	o	o	o	o			
Flex loop small	o	o	o	o	o	o	o			
Loop medium										
Endorectal	o	o	o	o	o	o	o			
TMJ right	o	o	o	o	o	o	o			
TMJ left	o	o	o	o	o	o	o			
CP peripheral angio	o	o	o	o	o	o	o			
CP breast array	o	o	o	o	o	o	o			
CP wrist										
Shoulder array	o	o	o	o	o	o	o			
CP extremity	o	o	o	o	o	o	o			
CP body array	o	o	o	o	o	o	o			
CP neck array	o	o	o	o	o	o	o			
CP spine array	o	o	o	o	o	o	o			
Head array	o	o	o	o	o	o	o			

SW Product No.:	024	087L	087S	047	089	085	086	007	083	091
Head CP (TX/RX)	-	-	-	-	-	-	-	-	x	x
Cardiac array	o	o	o	o	o	o	o			
MP 11								o		
MP 16								o		
MP 21s								o		
MP 28								o		
Head neck array								o		
Body spine array								o		

o = Option

8.17 Dicom

8.17.1 Dicom and UNIX nodes

For UNIX-based SIEMENS systems, such as NUMARIS 3x nodes, MagicView 1000 etc., make sure to define separate logical names for each service, e.g. Storage or Query & Retrieve.

8.17.2 Dicom configuration data for SIEMENS products

	Product	DICOM Service	AET	port no.
CT	Somaris 4: Somatom Plus4 Somatom AR (Spiral) Somaris 4.5: Balance	Storage SCU (send)	<PLA prefix> + DFOS	-
		Storage SCP (receive)	<PLA prefix> + DIC1	50082
		Query/Retrieve SCU	<PLA prefix> + DQRU	-
		Query/Retrieve SCP	<PLA prefix> + DQRY	50089
		Basic Print SCU	<PLA prefix> + DBPR	-
		Modality Worklist SCU	<PLA prefix> + DBWM	-
CT	Somaris 5: Somatom VZ	Storage SCU (send)	STU_<hostname> ¹	-
		Storage SCP (receive)	STU_<hostname>	104
		Query/Retrieve SCU	STU_<hostname>	-
		Query/Retrieve SCP	STU_<hostname>	104
		Basic Print SCU	MERGE_PRINT_SCU	-
		Modality Worklist SCU	HRI_<hostname>	-
CT	3D Virtuoso VA30	Storage SCU (send)	<HOSTNAME> ²	-
		Storage SCP (receive)	<HOSTNAME>	4444
		Query/Retrieve SCU	<HOSTNAME>	-
		Query/Retrieve SCP	<HOSTNAME>	4444
		Basic Print		

1. These are the default values for the AET. It is recommended to use the <hostname> for the AET. Future software versions will change the default value to the <hostname> as well.
2. Host name in uppercase letters

	Product	DICOM Service	AET	port no.
MR	Numaris 3: Magnetom Impact Magnetom Expert Magnetom Vision Magnetom Open Numaris 3.5: Harmony Symphony	Storage SCU (send)	<PLA prefix> + DFOS	-
		Storage SCP (receive)	<PLA prefix> + DIC1	50082
		Query/Retrieve SCP	<PLA prefix> + DQRY	50089
		Basic Print SCU	<PLA prefix> + DBPR	-
		Modality Worklist SCU	<PLA prefix> + DBWM	-
MR	syngo MR (Numaris 4)	Storage SCU (send)	STU_<hostname> ¹	-
		Storage SCP (receive)	STU_<hostname>	104
		Query/Retrieve SCU	STU_<hostname>	-
		Query/Retrieve SCP	STU_<hostname>	104
		Basic Print SCU	MERGE_PRINT_SCU	-
		Modality Worklist SCU	HRI_<hostname>	-

1. These are the default values for the AET. It is recommended to use the <hostname> for the AET. Future software versions will change the default value to the <hostname> as well.

	Product	DICOM Service	AET	port no.
SHS	MagicView 1000 VA31 / VB32	Storage SCU (send)	<PLA prefix> + DRSP ¹	-
		Storage SCP (receive)	<PLA prefix> + DRSP	50082
		Query/Retrieve SCU	<PLA prefix> + DQRY	-
		Query/Retrieve SCP	<PLA prefix> + DQRY	50089
		Basic Print SCU	<PLA prefix> + DBPR	-
SHS	MagicStore	Storage SCU (send)	<PLA prefix> + DRSP	-
		Storage SCP (receive)	<PLA prefix> + DRSP	50082
		Query/Retrieve SCP	<PLA prefix> + DQRY	50089 ²
		Query/Retrieve SCP	<PLA prefix> + DPQR	50090 ³
		Basic Print SCU	<PLA prefix> + DBPR	-

	Product	DICOM Service	AET	port no.
SHS	MagicView 50 / 300	Storage SCU (send) Storage SCP (receive) Query/Retrieve SCU	configurable ⁴ configurable configurable	- 104 -
SHS	MagicLink D / MagicServe	Storage SCU (send) Storage SCP (receive)	<PLA prefix> + DFOS <PLA prefix> + DIC1-8 ⁵	- 50082

1. Configurable with login: install -edit -config
2. Short Term database (IMS or IM)
3. Long Term database (PDIR or MM)
4. it is recommended to use the <hostname> as AET
5. 8 import channels (DIC1 to DIC8) with different PACSnet destinations can be configured

Migration from Frame-Word to CMS.

Chapter 5: Information added how to Backup/Restore site specific data on USB device.

There are no Hazard IDs in this document.

All documents may only be used for rendering services on Siemens Healthcare Products. Any document in electronic form may be printed once. Copy and distribution of electronic documents and hardcopies is prohibited. Offenders will be liable for damages. All other rights are reserved.

healthcare.siemens.com/services

Siemens Healthineers Headquarters
Siemens Healthcare GmbH
Henkestr. 127
91052 Erlangen
Germany
Telephone: +49 9131 84-0
siemens.com/healthineers

Print No.: MR-000.816.17.05.02 | Replaces: MR-000.816.17.04.02
Doc. Gen. Date: 07.11 | Language: English
Part No.: n.a.
© Siemens Healthcare GmbH, 2004

siemens.com/healthineers