


TeamViewer Instructions


Download version 12.xx at <http://www.teamviewer.com>

Run File:

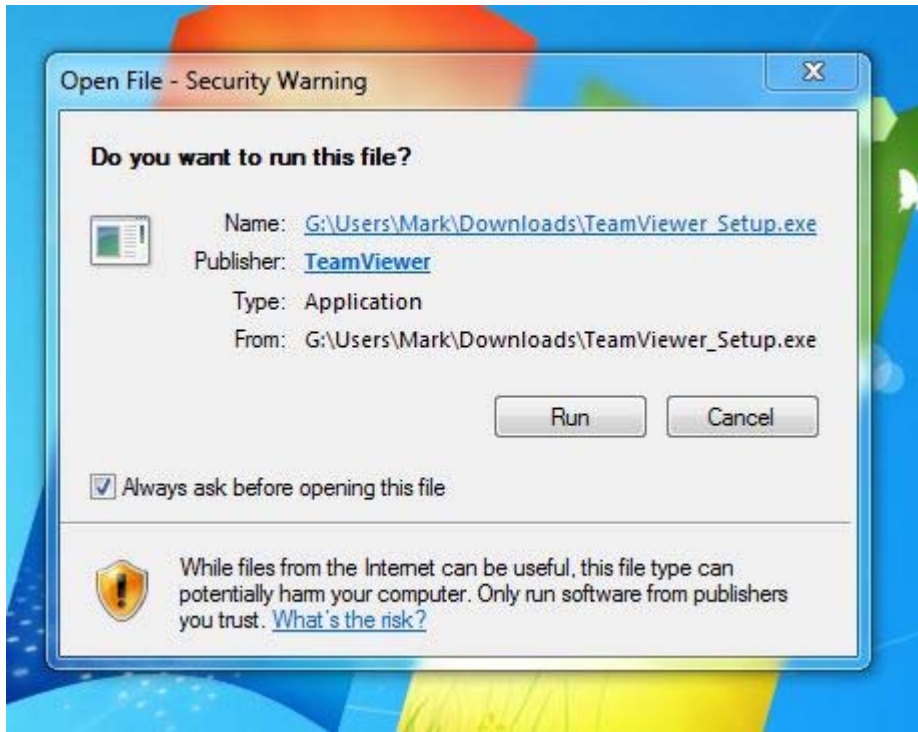

Check Basic Installation if you have Admin Rights on your computer, otherwise check Run Once

Please note you will have to run the exe file every time you need to use TeamViewer if you do not have Admin rights.

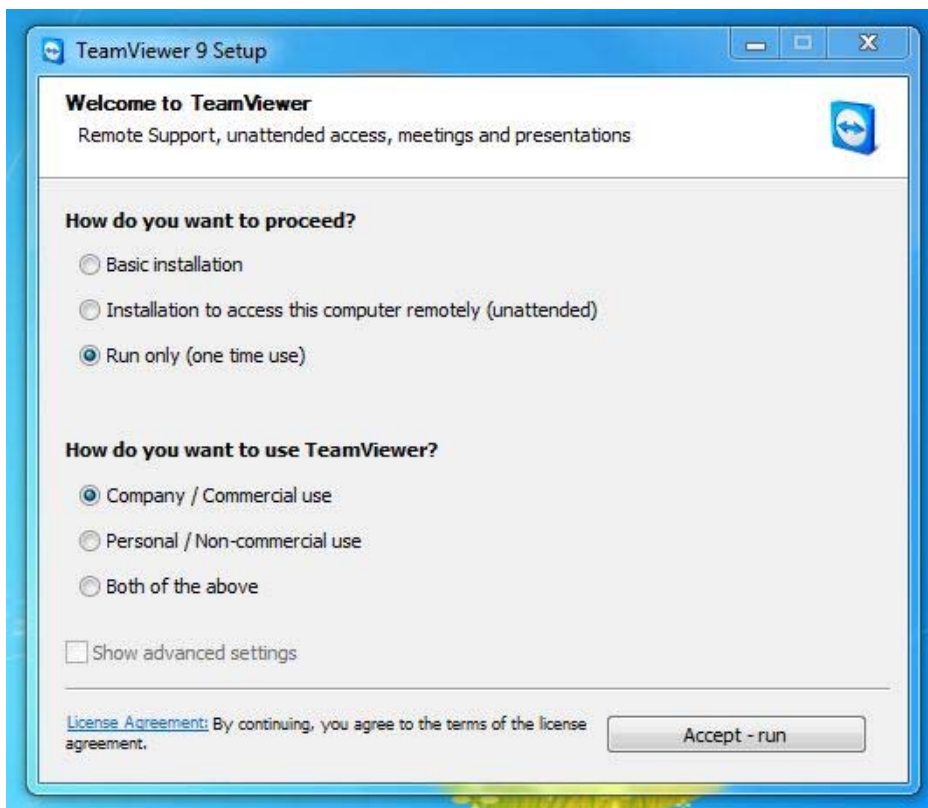

Enter email address and password then click Sign In

email: bct_engineer@coldtrac.net

Password: [cOldtrac\\$](#)

The image shows the TeamViewer application interface. The main window is titled 'TeamViewer' and has a menu bar with 'Connection', 'Extras', and 'Help'. Below the menu bar, there are two tabs: 'Remote Control' and 'Meeting'. The 'Remote Control' tab is active. The interface is divided into two main sections: 'Allow Remote Control' and 'Control Remote Computer'. The 'Allow Remote Control' section displays the user's ID as '301 319 558' and their password as '978xkg'. It also includes a 'Personal Password' field with a lock icon. The 'Control Remote Computer' section has a 'Partner ID' field and radio buttons for 'Remote control' (selected) and 'File transfer'. A 'Connect to partner' button is located at the bottom of this section. On the right side, there is a 'Computers & Contacts' window titled 'Please sign in...'. It contains an 'E-mail' field with the address 'bct_engineer@coldtrac.net', a 'Password' field with masked characters, and a 'Sign In' button. Below the 'Sign In' button, there is a checkbox for 'Keep me signed in' and two links: 'Forgot password' and 'Sign Up'. At the bottom of the TeamViewer window, a status bar indicates 'Ready to connect (secure connection)' and a 'Computers & Contacts' button with a double arrow icon.

TeamViewer

Connection Extras Help

Remote Control Meeting

Allow Remote Control

Please tell your partner the following ID and password if you would like to allow remote control.

Your ID 301 319 558

Password 978xkg

Use your personal password to access this computer from anywhere.

Personal Password

Control Remote Computer

Please enter your partner's ID in order to control the remote computer.

Partner ID

☒ Remote control
☐ File transfer

Connect to partner

Computers & Contacts

Computers & Contacts

Please sign in...

E-mail

bct_engineer@coldtrac.net

Password

Sign In

☐ Keep me signed in


[Forgot password](#)

[Sign Up](#)


Ready to connect (secure connection)

All ColdTRAC Sites will be shown in list divided by Organization.

Double Click on site you wish to view.

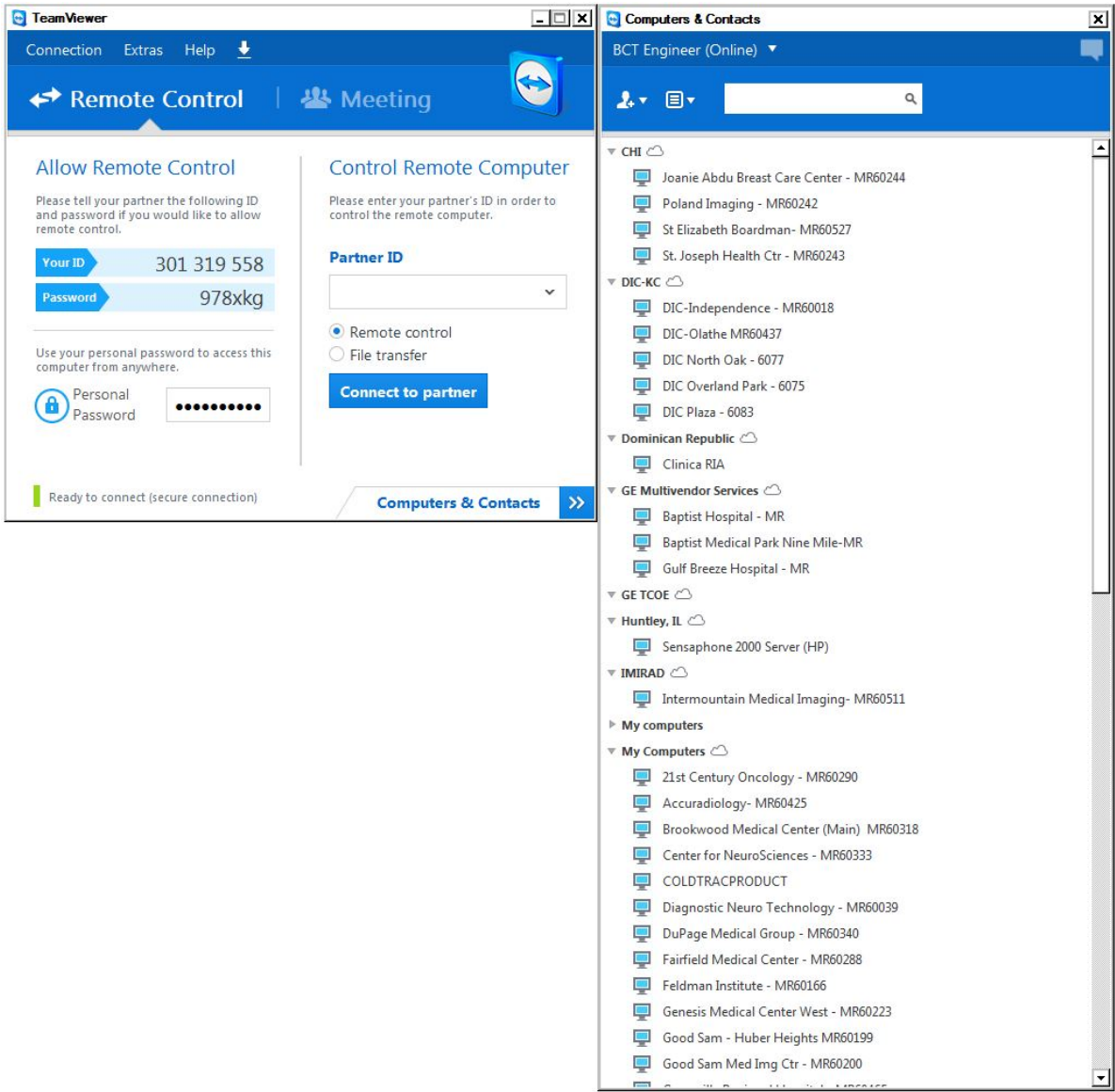

You will then see the ColdTRAC App screen.

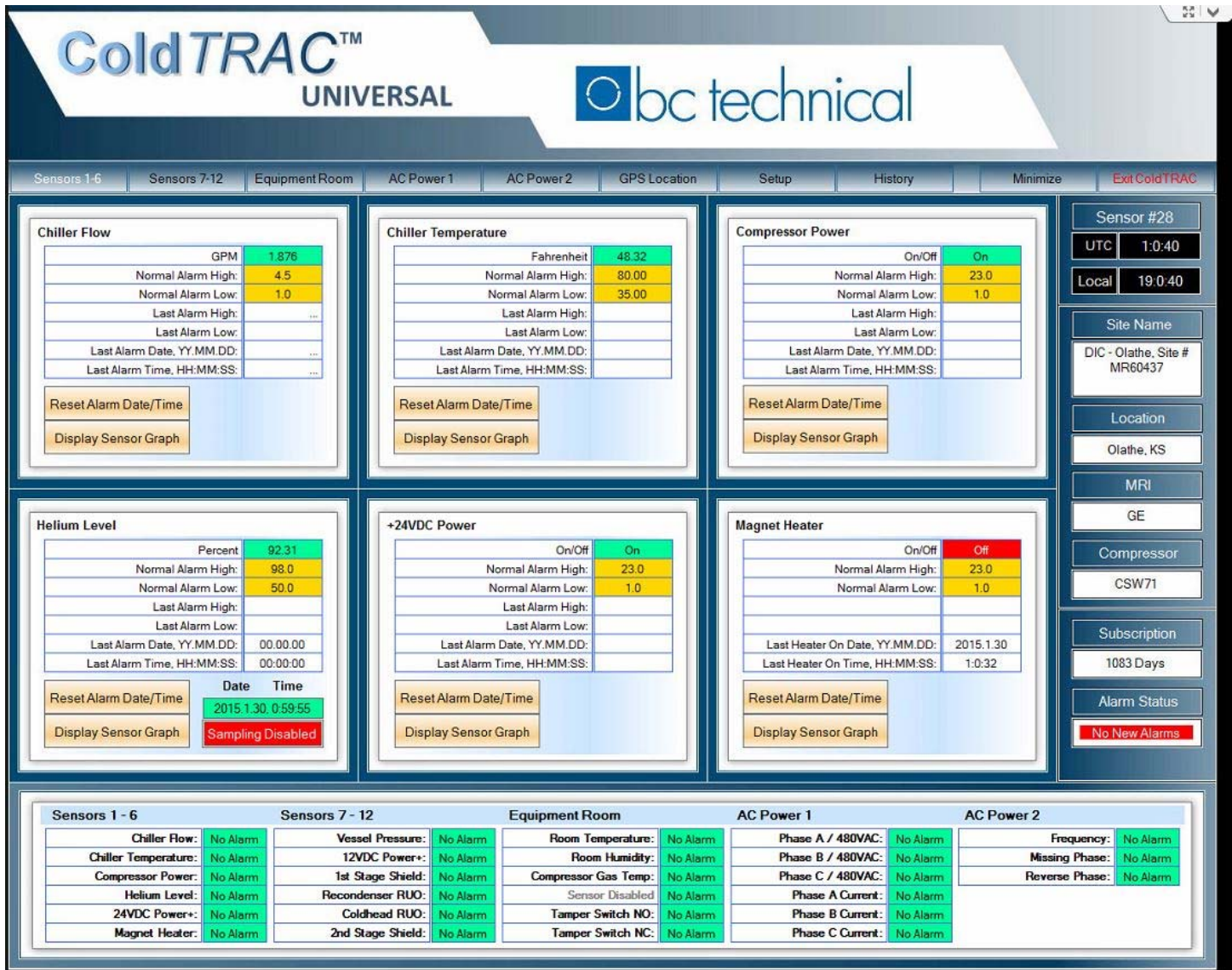

Click on the tabs across the top to view all sensor values. Please do not click on GPS location, Setup, or History as these are not end user viewable or for future upgrades to our software.

If you need to turn off or suspend an alarm please call or email one of the ColdTRAC Admin Team. We will suspend the alarm and set a reminder appointment in Outlook to reinstate alarm. The Engineer will be included in the appointment.

To exit the session click on the X in the Teamviewer Bar at the top of the screen, or close the remote browser window.

The screenshot shows the ColdTRAC UNIVERSAL interface. At the top, there is a navigation bar with the ColdTRAC UNIVERSAL logo and the bctechncal logo. Below the navigation bar, there are several tabs: Sensors 1-6, Sensors 7-12, Equipment Room, AC Power 1, AC Power 2, GPS Location, Setup, History, Minimize, and Exit ColdTRAC. The main content area is divided into three columns, each displaying a different set of data:

- Chiller Flow:** A table showing GPM (2.808), Normal Alarm High (4.5), Normal Alarm Low (1.0), Last Alarm High (0.407), Last Alarm Date (2015.8.27), and Last Alarm Time (18:26:18). Below the table are buttons for "Reset Alarm Date/Time" and "Display Sensor Graph".
- Chiller Temperature:** A table showing Fahrenheit (64.28), Normal Alarm High (80.00), Normal Alarm Low (35.00), Last Alarm High (81.15), Last Alarm Date (2015.8.27), and Last Alarm Time (18:27:5). Below the table are buttons for "Reset Alarm Date/Time" and "Display Sensor Graph".
- Compressor Power:** A table showing On/Off (On), Normal Alarm High (23.0), Normal Alarm Low (1.0), Last Alarm High (Off), Last Alarm Date (2015.9.2), and Last Alarm Time (0:54:59). Below the table are buttons for "Reset Alarm Date/Time" and "Display Sensor Graph".

On the right side of the interface, there is a sidebar with a "Sensor #6" section showing UTC (1:38:19) and Local (20:38:19) time. Below this is a "Site Name" section showing "St Elizabeth Boardman Health Center # MR60527". At the bottom of the sidebar is a "Location" section showing "Boardman, OH" and "MRI".

Once closed the audit window will appear. Please enter what you were looking at and your name and hit save.

The screenshot shows a "Comment" window. The window has a title bar with a "Comment" label and a close button. The main content area has a blue header with the word "Comment". Below the header, there is a text area with the following text:

Please comment on your connection with St Elizabeth Boardman- MR60527 at 09/01/2015 8:39 PM

Below this text is a text input field containing the text "Looked at chiller values" and "M Abbs". At the bottom of the window, there is a "Nothing to comment" link and a "Save" button.

If at any time you have questions or problems logging in or viewing the App please call or email

ColdTRAC Admin Team:

Mark Abbs 913-489-1320

mark.abbs@bctechncal.com

Rob Mueller 847-533-1603

rob.mueller@bctechncal.com