

MAGNETOM

Update Instructions Evolve & Enhancement - MR002/18/E

VE11D Evolve

Update within deadline	12 months		
Mandatory	☐ Yes	☒ No	
Affected	☐ Warranty	☒ Contract	☐ Others
Update by	☐ Customer	☐ Informatics	☐ RSC ☐ Apps ☒ CSE
Remote update	☐ Yes	☒ No	
Intranet download available	☐ Yes	☒ No	
Update material required	☒ Yes	☐ No	
Material free of charge	☒ Yes	☐ No	☐ No, credit if returned
Return of parts	☐ Yes	☒ No	
Application training needed	☒ Yes	☐ No	☐ Recommended
Application training time	see Customer application training duration		
Estimated completion time	7h	Number of CSEs: 1	

Remarks: Together with the Evolve installation date the customer application training has to be coordinated!

Scope

Material number	10684331, 10684332, 10684333, 10684334
Software version	VD13A
Related to Update Instructions	MR015/18/E VE11D Application Training

Change reference no.: 672407
 Name: D. Hofmann
 Department: SHS CS SLM MR SW&OP

Document Version

Siemens Healthcare GmbH reserves the right to change its products and services at any time.

In addition, manuals are subject to change without notice. The hardcopy documents correspond to the version at the time of system delivery and/or printout. Versions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens Healthineers office to order a current version or refer to our website <http://www.siemens-healthineers.com>.

Disclaimer

Siemens Healthcare GmbH provides this documentation “as is” without the assumption of any liability under any theory of law.

The content described herein shall be used by qualified personnel who are employed by Siemens Healthcare GmbH or one of its affiliates or who are otherwise authorized by Siemens Healthcare GmbH or its affiliates to use such documents.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens Healthcare GmbH or one of its affiliates are not entitled to use this documentation without prior written authority.

Copyright

“© Siemens Healthineers, 2018” refers to the copyright of a Siemens Healthineers entity such as:

Siemens Healthcare GmbH - Germany

Siemens Shenzhen Magnetic Resonance Ltd. - China

Siemens Shanghai Medical Equipment Ltd. - China

Siemens Healthcare Private Ltd. - India

Siemens Medical Solutions USA Inc. - USA

Siemens Healthcare Diagnostics Inc. - USA and/or

Siemens Healthcare Diagnostics Products GmbH - Germany

1	Preparation	5
1.1	General Information	5
1.1.1	Systems/Products Affected	5
1.1.2	Reason for the Update	5
1.1.3	syngo Feature	5
1.1.4	Customer application training duration	6
1.1.5	Prerequisites	6
1.1.6	Special Tools	7
1.2	Material Information	8
1.2.1	Ordering Information	8
1.2.2	Content of the Update Kit	8
1.2.3	Return of Parts	9
2	Evolve VE11D - First Steps	10
2.1	First work steps for the VE11D Evolve upgrade	10
2.1.1	Overview	10
2.1.2	License-relevant Data	11
2.1.3	MRWP (MRSC) Step: Backing up Dynamic Data at the MRWP (MRSC)	12
2.1.4	MRAWP (MRC) Step: Backing up Dynamic Data at the MRAWP (MRC)	14
2.1.5	Sending the Backup CDs to CV ME	18
2.1.6	Perform Remote Backup for the Evolve Update	19
3	Evolve VE11D - Final Steps	37
3.1	Final Work Steps for EVOLVE from VD13A to VE11D	37
3.1.1	Overview	37
3.1.2	MRWP (MRSC) Preparatory Steps for EVOLVE Update at the MRWP with the Old Software Version	38
3.1.3	MRAWP (MRC) Preparatory Steps for EVOLVE Update at the MRAWP with the Old Software Version	39
3.1.4	Replacing the MRAWP/MRWP and MARS	41
3.1.5	MRAWP (MRC) System Start-up	46
3.1.6	MRAWP (MRC) New secure Service Key Upload	47
3.1.7	MRAWP (MRC) Configuration checks	48
3.1.8	MRAWP (MRC) Tune-up and QA Steps	53
3.1.9	MRAWP (MRC) Surface Coils QA Steps	55
3.1.10	MRAWP (MRC) Network Printer Installation	55
3.1.11	MRAWP (MRC) Function Test of the DICOM Nodes	55
3.1.12	MRAWP (MRC) System Management	55
3.1.13	MRAWP (MRC) Restoring Changed Customer Measurement Programs from VD13A	61
3.1.14	MRAWP (MRC) Manual Restore of Voice Commands	63
3.1.15	MRAWP (MRC) Standby Function Check	64
3.1.16	MRAWP (MRC) Setup Advanced User Login Function	65
3.1.17	MRAWP (MRC) Setup Auto Login Function	65
3.1.18	MRAWP (MRC) Performing COBR Image backup	65
3.1.19	MRAWP (MRC) Enable the Auto Login Function	66

3.1.20	MRWP (MRSC) Startup	67
3.1.21	MRWP (MRSC) Checking the Patient Browser Connection with Shared Database	68
3.1.22	MRWP (MRSC) Network Printer Installation on the MRSC	68
3.1.23	MRWP (MRSC) Configuration checks	69
3.1.24	MRWP (MRSC) Function Test of the DICOM Nodes	73
3.1.25	MRWP (MRSC) System Management	73
3.1.26	MRWP (MRSC) Setup Auto Login Function	79
3.1.27	MRWP (MRSC) Setup Advanced User Login Function	80
3.1.28	MRWP (MRSC) Perform a COBR Image Backup on the MRWP (MRSC)	80
3.1.29	MRWP (MRSC) Enable the Auto Login Function	80
3.1.30	Destroying Data on Old Host PCs	81
3.1.31	Final Check - to be performed for MRAWP (MRC) and MRWP (MRSC) (if applicable)	84
3.1.32	Customer Information	87
3.1.33	Removed Windows Feature in VE11D	87
3.1.34	Final Work Steps	87
3.1.35	Changes to Previous Version	88
4	Additional Information	89
4.1	Appendix	89
4.1.1	Service key handling with new Secure Service Key	89
4.1.2	New Password Handling in syngo	94
4.1.3	Advanced User handling	98
4.1.4	New feature for DICOM node configuration	100
4.1.5	Drive Letter Assignment for the CD/DVD-ROM Drive and/or for the CD/DVD-R Drive	101
4.1.6	Correct syngoUI Language Settings	103
4.1.7	Failure messages where the regional language setting did not match, e.g. UI language	105
4.1.8	Language settings	107
4.1.9	System management installation workaround	139
4.1.10	Set Embedded Control in update mode	146
4.1.11	Checking and setting the MOXA - Ports for HPz440	148
5	Completion Protocol MR002/18/E	159
5.1	VE11D Evolve	159
5.1.1	Reason for the Update	159
5.1.2	Protocol	160
6	List of Hazard IDs	161

1.1 General Information

1.1.1 Systems/Products Affected

See page 1 and a valid Evolve contract has to be assigned to the system.

1.1.2 Reason for the Update

This Update Instruction is for obligating Evolve update. It provides the compatible hardware that is required (host(s) and MaRS) and the software syngo MR E11D according to the customer's Evolve agreement. The software VE11D does contain:

- Significant enhancement in security configuration for secure scanner operation in hospital network.
- syngo Security Enhanced

In case the previous SW version was VD13 the VE11D provides additionally the new Dot Cockpit for easy Dot Engine configuration together with several workflow and performance enhancements:

- DotGo
- Dot Cockpit

With the available options:

- StarVibe (FREEZEit)
- TWIST-Vibe (FREEZEit)
- LiverLab
- MyoMaps
- Dot Engines (Spine, Large Joint, Breast)

1.1.3 syngo Feature

VH22D

- With version VH22A, syngo introduced several measures recommended by CERT and DoD to improve system IT security. Password handling has changed completely, which could impact field service significantly. There are no longer fixed or default passwords within syngo. Only random or individually defined passwords are valid. For new password handling, see (→ New Password Handling in syngo / Page 94) and for the Advanced User from MR, see (→ Advanced User handling / Page 98).

syngo Expert-i

- Interactive remote assistance during the MR exam. Get real-time access to imaging data and exam information from any networked PC within the hospital network. Until

now, radiologists or other experts had to stop what they were doing and go to the MR scanner when complex clinical questions arose.

Now, questions can be addressed quickly and efficiently via a networked PC.

Some of the benefits of syngo Expert-i are:

- » Excellent results right from the first examination
- » Streamlined workflow & faster patient throughput
- » Reduces repeat rates with a check on images while the patient is still in the examination room
- » Reduces training effort by enabling expert assistance for specialized procedures

1.1.4 Customer application training duration

■ Upgrade for MAGNETOM Avanto and Verio from MR D13 to MR E11:

- » **5 days**
In case of Upgrade and the staff is not familiar with E11 software
- » **3-5days**
In case the staff is already familiar with MR E11 software



Clarify the customer situation to plan the appropriate on-site training duration.

1.1.5 Prerequisites

The latest version of this document can be found in the Services Knowledge Base (SKB).



Be aware that with VE11D, general syngo security features are activated and that the account will be locked if the password is entered incorrectly four times. Please complete the syngo Security WBT.



"Work In Progress" packages provided for a previous software version (VD13*) will generally NOT work with VE11D!



Due to DICOM changes effective with VE11D, adaptations regarding the locally connected DICOM nodes may become necessary on site.



The BreVis evaluation SW is already included in VE11D. It is no longer necessary to order a separate installation medium.

- To also use IDEA with MR systems running 4VE11D, a current version of the IDEA software is necessary.
 - Request the current IDEA package from HC DI MR PI CM for sites where this feature is being used.
 - Please note: In accordance with Siemens policy, IDEA will be provided within the framework of a cooperation agreement only. For questions, please contact MR co-operation management - HC IM MR PI CM.
- If a Leonardo/MMWP workstation (MultiModality Workplace) is used to process data acquired with the MR system with NUMARIS/4 E11, the current Leonardo/MMWP workstation (MultiModality Workplace) versions are not compatible. The new Leonardo/MMWP workstation (MultiModality Workplace) version MM WP VE61A is required.

The following host PC types are supported by MMWP version VE61A:

 - R670-2
 - M470-2
 - M720

If the customer has a different host PC, the hardware must be exchanged as well.
- Ensure that the VE11D application training timing is suitable for your site.
- In case the customer uses 3rd party coils on the MR System which require the installation of coil files, please contact the coil manufacturer directly or ask the customer for it.
- With VE11D the new Secure Service Key will be provided, for more details please refer to (→ Service key handling with new Secure Service Key / Page 89).

1.1.6 Special Tools

n.a.

1.2 Material Information

1.2.1 Ordering Information

The following update kit has to be ordered via the METO (Medical Electronics Tracking and Ordering) process.

Place the order

- directly in METO at Siemens Healthcare / CV and
- at Siemens Healthcare / CS ML (SAP Distribution Channel, Plant 2050)

To order the update via the web application, the following information has to be provided:

- System type
- System serial number
- Keyboard layout (German, US, French, Spanish, Swedish,...)
- Shipping address for the update
- Desired delivery date for the upgrade
- For countries with more than one language, the preferred language of the customer

1.2.2 Content of the Update Kit

The material delivered for a specific system depends on the configuration. Below you can find the most important materials that, if applicable, will be delivered with an upgrade:

Material No.	Name	MRAWP only	MRAWP + MRWP
		UI Kit# 11344763	UI Kit# 11344764
11061964	MR HWCfg Container	1	1
10432853	DVD COBR-Recovery-Image	1	2
11344758	Expert-I	1	1
07752103	Set of blank DVDs	1	2
	Customer documentation package	1	1
	MaRS	1	1
	Keyboard	1	2
	Mouse	1	2
	Host (MRAWP)	1	1
	Host (MRWP)	0	1
	Cable	1	2

1.2.3 Return of Parts

Please return the hardware that was replaced by Evolve if the customer agrees that Siemens should dispose of the old HW.

Before proceeding with the disposal, please make a note on the last page of the UI (Completion Protocol) e.g.: Siemens will handle the disposal of the material with part number: ...(list with material and serial number for disposal)...

This note has to be signed by the customer stating that he agrees to Siemens' handling the disposal of the listed hardware.

- All returns (including when receiving defective or unused materials or repair parts) are to be handled in accordance with your ARTD **"Principles of Return Goods Management"**.
- Returns have to be properly packed according to local **"Transport, Packaging and Storage of Spare Parts"** regulations.
- Use the same boxes that were used to pack the Evolve hardware
- Do not change the Evolve labels on the boxes. This way CS ML can identify the return as an Evolve return (if necessary, write the Evolve update material number on the boxes).


For disposal or if the customer wishes to keep the hardware, you must perform the work steps in the chapter: "Destroying Data on Old Host PCs"

2.1 First work steps for the VE11D Evolve upgrade

2.1.1 Overview


2.1.1.1 Performing an On-site Backup for the Evolve Upgrade

Fig. 1: Flowchart - local backup

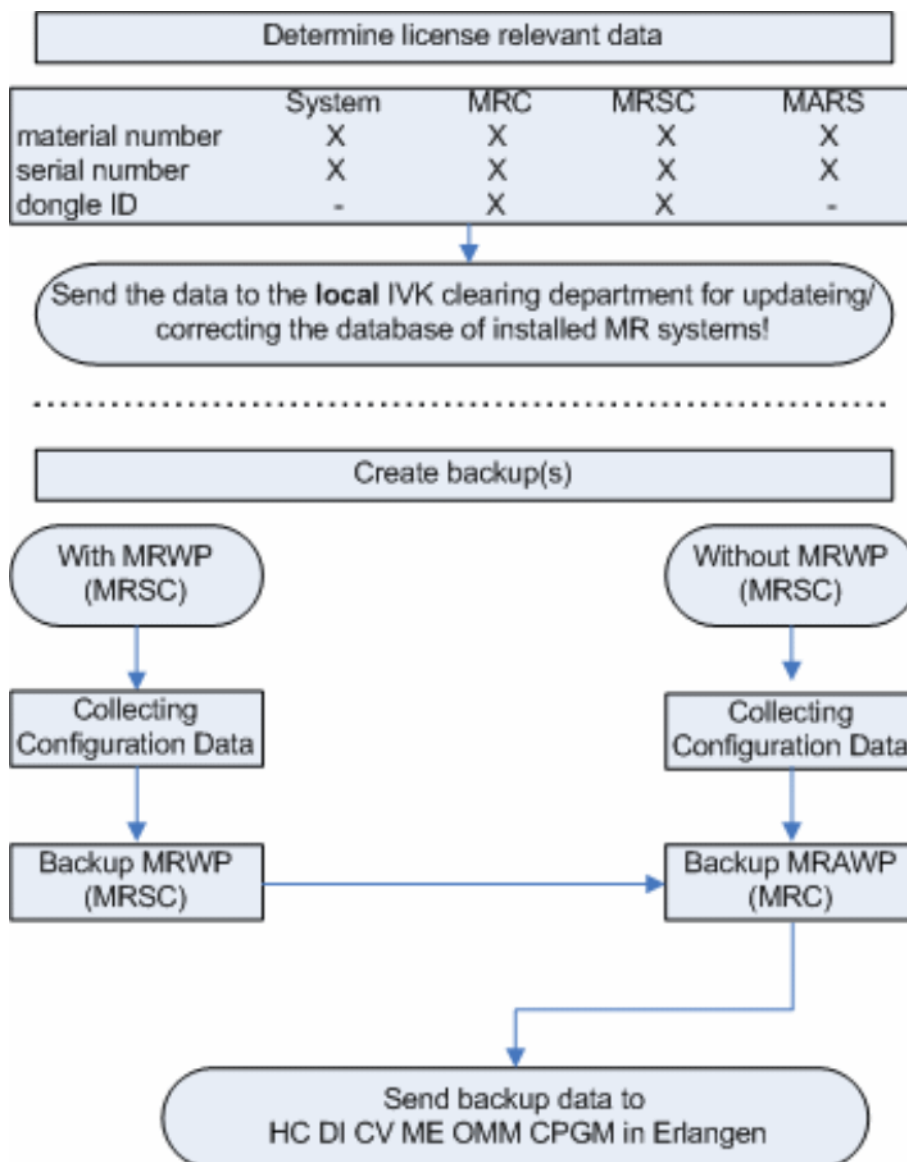

2.1.1.2 Performing a Remote Backup for the Evolve Upgrade

Fig. 2: Flowchart - remote backup

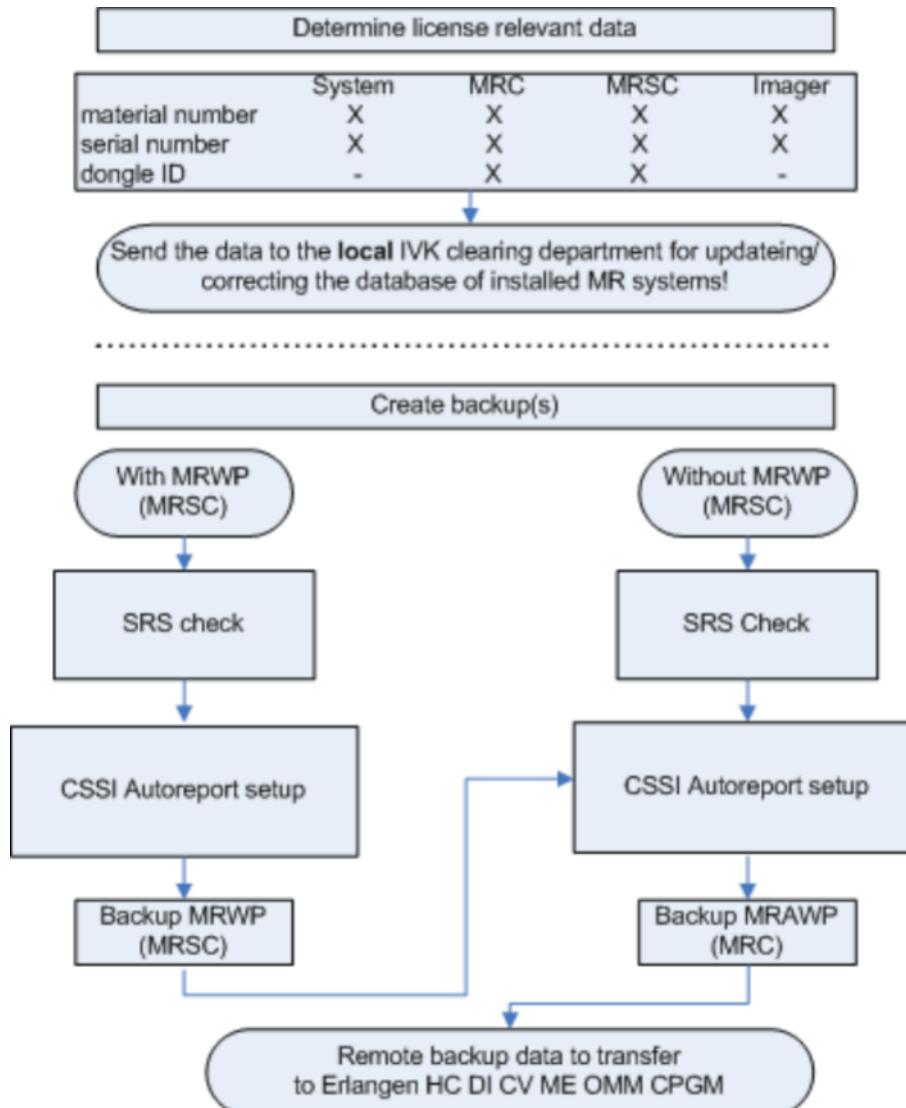

2.1.2 License-relevant Data

The new licensing application (Licensing@Healthcare) is dependent on accurate IVK data populated via the local databases of the RUs.

The current data quality is inadequate and in most cases does not allow valid license files to be created for CSSI- and Evolve upgrades.

For this reason, the following data must be determined on site and forwarded to the IVK clearing of the RU for checking/correcting the relevant data in the local database before a CSSI or Evolve upgrade is ordered.

MRAWP (MRC)	Material#		From the stickers on the components showing "IVK".
	Serial#		
MRWP (MRSC) (if applicable)	Material#		
	Serial#		
Imager	Material#		
	Serial#		

The dongle IDs which are also important are listed in the license.dat files of AWP and MRWP, which are part of the corresponding backups to be sent to Erlangen for any CSSI or Evolve upgrade.

2.1.3 MRWP (MRSC) Step: Backing up Dynamic Data at the MRWP (MRSC)



If you have a satellite console with a shared database connected to the main console, perform the following steps first.

If no MRWP (MRSC) is present, skip the following steps and start at (→ MRAWP (MRC) Step: Backing up Dynamic Data at the MRAWP (MRC) / Page 14)

2.1.3.1 Collecting MRWP Configuration Data from Local Service

ACTION	MRWP
■ Select Options --> Service --> Local Service .	☐
■ Enter the valid service key and click OK .	☐
■ Select Reports in the service platform.	☐
■ On the Report Files page, select Site Info Imprint	☐
■ Print out the Imprint report (generate a PDF file or write down the system options) Launch the update of the Imprint file (bottom right corner) In order to generate a PDF file of the Imprint report, press " show as PDF ". Save the PDF file under C:\MedCom\MriSiteData\Imprint.pdf	☐
■ On the Report Files page, select Site Info Configuration	☐

ACTION	MRWP
- Print out the Configuration report (generate a PDF file or write down the system options) - DICOM networking (write down all DICOM network nodes and settings) - DICOM print devices (write down all DICOM print devices and settings) - DICOM HIS/RIS (write down all DICOM HIS/RIS nodes and settings) In order to generate a PDF file of the Configuration report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\Report.pdf	☐
On the Report Files page, select Site Info MR Configuration	☐
Print out the MR Configuration report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\MrReport.pdf	☐
Click EXIT to close the Report Files popup window	☐



It is recommended to copy the PDF Report Files to an external medium to have the information for the configuration check in the new software version.

Please check the external medium first with the virus scanner of the service PC. Ensure that the external medium used is free of viruses.

2.1.3.2 Master Backup at the MRWP (MRSC)

ACTION	MRWP
Perform a Master Backup in local service. The Master Backup includes the Numaris4 Package, Software Settings, and Security Settings. Label an empty CD-R medium with Systemserial# Master Backup VD13A MRWP (actual date) Start local service, select Backup&Restore, insert the labeled CD-R medium into the [CD-R] drive (S:\), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu, select Master-Backup in the Package menu, start the backup by clicking OK.	☐

ACTION	MRWP
Note: - Check the output log for errors. If any other error occurs apart from those mentioned below, the master backup has to be created again. - In some cases, the following popup message: "Please insert the media into Drive S:\ and press 'OK' to proceed " appears. Click "OK" to continue - If the following error messages are displayed at the end of the backup routine, ignore them. The backup is then completed correctly. - Error Execute starting components... 3 rti.exe runComp an_server rti->runComp(an_server, Pid, Container) failed and/or - Error Execute starting components... 2 rti.exe runComp stu_server rti->runComp(stu_server, Pid, Container) failed and/or - Error Execute starting components... 1 rti.exe runComp ima_server rti->runComp(ima_server, Pid, Container) failed ...	
Remove the CD.	☐
Click Home and close the local service UI.	☐



If the backup creation to CD fails, please create the master backup to hard disk and transfer it to your service PC. On your service PC, burn the master backup to CD.

2.1.4 MRAWP (MRC) Step: Backing up Dynamic Data at the MRAWP (MRC)





If you have a satellite console with a shared database connected to the main console, start with the action item:

- (→ MRWP (MRSC) Step: Backing up Dynamic Data at the MRWP (MRSC) / Page 12)

This item has to be performed first. If this step has already been completed on the MRWP, continue with the MRAWP.

2.1.4.1 Collecting MRAWP Configuration Data from Local Service

ACTION	MRAWP
Select Options --> Service --> Local Service.	☐
Enter the valid service key and click OK.	☐
Select Reports in the service platform.	☐

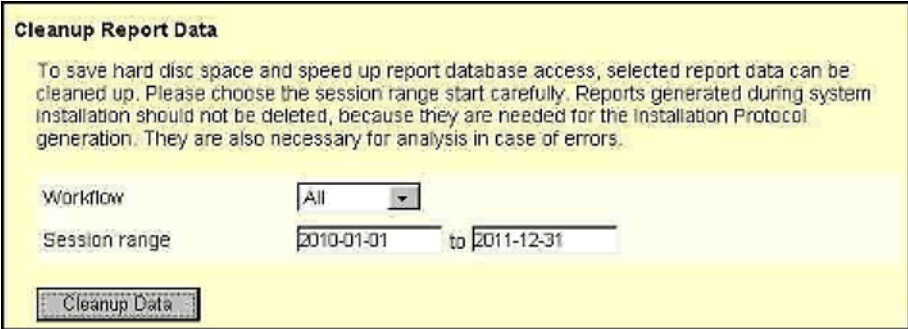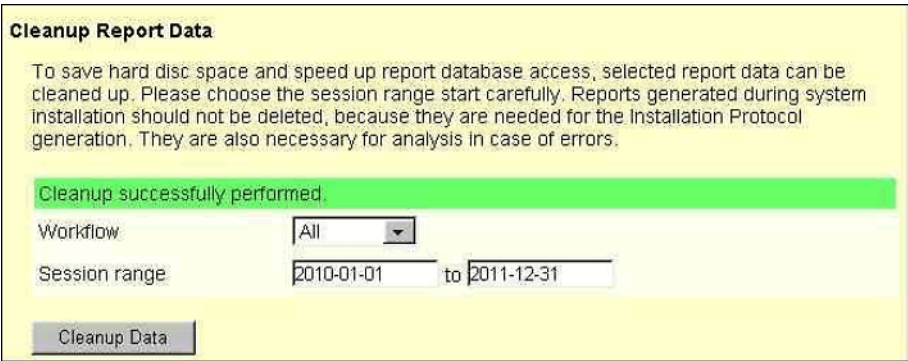ACTION	MRAWP
■ Purge service-related files (test tool files and tune-up, QA files older than one year). If the files are not purged, the size of the master backup created will be huge. - Select Reports -> Clean up - Under Cleanup Report Data select Workflow All - Session range 2010-01-01 to (today - 1 year) e.g. 2011-07-01 - You can also click in the fields and a calendar popup appears. <i>Fig. 3: Data Cleanup</i>  (1) Start date (2) End date - Press the Cleanup Data button. <i>Fig. 4: Confirmation popup</i>  Confirm the popup with OK. - The following shows a successful file purge. <i>Fig. 5: Successful cleanup message</i> 	<input data-bbox="1385 275 1425 315" type="checkbox"/>
■ On the Report Files page, select Site Info Imprint	<input data-bbox="1385 1821 1425 1861" type="checkbox"/>


ACTION	MRAWP
Print out the Imprint report (generate a PDF file or write down the system options) Launch the update of the Imprint file (bottom right corner) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\Imprint.pdf	☐
On the Report Files page, select Site Info Configuration	☐
- Print out the Configuration report (generate a PDF file or write down the system options) - DICOM networking (write down all DICOM network nodes and settings) - DICOM print devices (write down all DICOM print devices and settings) - DICOM HIS/RIS (write down all DICOM HIS/RIS nodes and settings) In order to generate a PDF file of the Configuration report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\Report.pdf	☐
On the Report Files page, select Site Info MR Configuration	☐
Print out the MR Configuration Report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\MrReport.pdf.	☐
On the Report Files page, select Site Info CAN24 Configuration	☐
Print out the CAN24 Configuration report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\CAN24Report.pdf	☐
On the Report Files page, select Site Info MBox Configuration	☐
Print out the MBox Configuration report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\MBoxReport.pdf	☐
Click EXIT to close the Report Files popup window	☐



It is recommended to copy the PDF report files to an external medium to have the information for the configuration check in the new software version.

Please check the external medium first with the virus scanner of the service PC. Ensure that the external medium used is free of viruses.

2.1.4.2 Master Backup at the MRAWP (MRC)

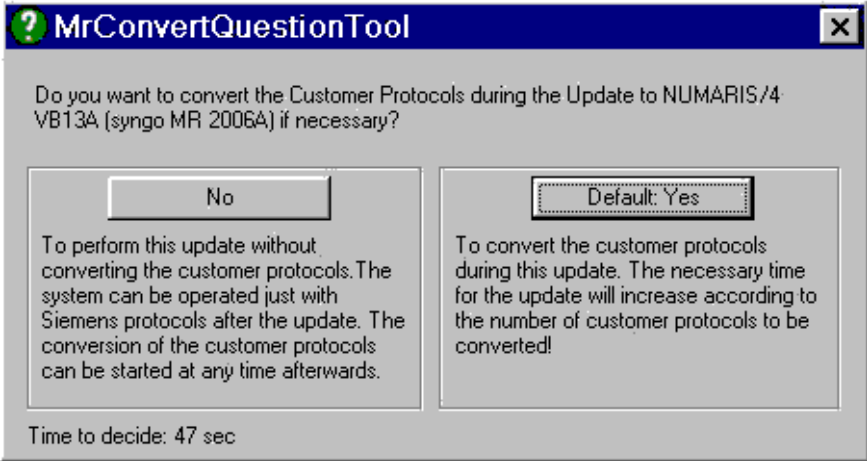ACTION	MRAWP
Perform a Master Backup in local service. The Master Backup includes the Numaris4 Package, Software Settings, and Security Settings. Label an empty CD-R medium with Master Backup VD13A MRAWP (actual date) Start local service, select Backup&Restore, insert the labeled CD-R medium into the [CD-R] drive (S:\), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu, select Master-Backup in the Package menu, start the backup by clicking OK. During the Master Backup, the following popup appears: Fig. 6:  Press the "Yes" button. At this point, this popup will not start the protocol conversion if "yes" was selected. This popup determines the default focus ("YES" or "NO") of the MrConvertQuestion-Tool when restoring the Master Backup.	<input data-bbox="1385 584 1422 622" type="checkbox"/>

ACTION	MRAWP
Note: ■ Check the output log for errors. If any other error occurs apart from those mentioned below, the master backup has to be created again. ■ In some cases, the following popup message: "Please insert the media into Drive S:\ and press 'OK' to proceed" appears. Click "OK" to continue ■ If the following error messages are displayed at the end of the backup routine, ignore them. The backup is then completed correctly. - Error Execute starting components... 3 rti.exe runComp an_server rti->runComp(an_server, Pid, Container) failed and/or - Error Execute starting components... 2 rti.exe runComp stu_server rti->runComp(stu_server, Pid, Container) failed .. and/or - Error Execute starting components... 1 rti.exe runComp ima_server rti->runComp(ima_server, Pid, Container) failed ...	
■ Remove the CD.	☐
■ Click Home and close the local service UI.	☐



If the backup creation to CD fails, please create the master backup to hard disk and transfer it to your service PC. On your service PC, burn the master backup to CD.

2.1.5 Sending the Backup CDs to CV ME

Send the backup to CV ME via electronic data transfer (Siemens Secure File Exchange (SecuFEx), Big File Transfer, internal transfer folder,...). Please always add the serial number of the system into the file name.

Evolve mail account: **immrevolve.team@siemens-healthineers.com**

If you send the backup via SecuFEx to our general Evolve mail account, please consider that the PKI cannot be applied to this general mail account, therefore please select "Low Security" option with a temporary password to provide access to the backup data!

In case with SecuFEx (4 steps):

E-Mail: **immrevolve.team@siemens-healthineers.com**

in step 3 of 4 select - Access Type: **Low Security** (access via temporary password)

Mail Subject: **EVOLVE/CSSI Customer Data**

CV contact:

Siemens Healthcare GmbH

SHS DI CV ME OMM CPGM

Attn.: Mrs. Sina Aue, sina.aue@siemens-healthineers.com

Allee am Roethelheimpark 2

D - 91052 Erlangen

2.1.6 Perform Remote Backup for the Evolve Update



Before the remote backup can be performed, you must inform the customer about the duration of the backup and the requirement for full access. Also ask the customer to keep the system running in standby mode for the backup transfer.

2.1.6.1 Prerequisites

- System is connected to the SRS infrastructure via a VPN connection
- System is in **full remote access mode**
- The estimated time is about 30-40 min. for the remote backup, depending on the amount of backup data.

2.1.6.2 SRS Configuration Check

- Remote login to the system using the cRSP (common Remote Service Platform)
 - Check the RDIAG function on the **MRAWP** and, if applicable, on the **MRWP** as well
 - Test the FTP connection to the RDIAG server as follows:
 - Local Service -> File&Image Tools
 - Source: File
 - Action: Transfer to System
 - File Type: Any file
 - From source: SRS_ftp
- If a list from the directory /outgoing/mr (directory on the RDIAG server) is displayed, the connection and FTP access are OK



If the FTP test fails, check the SRS configuration; refer to (→ Final Check - to be performed for MRAWP (MRC) and MRWP (MRSC) (if applicable) / Page 84)



If you do not have a satellite console with a shared database connected to the main console, the step (→ Auto Report Setup and Remote Backup on MRWP / Page 20) can be skipped.

In this case, continue with (→ Auto Report Setup and Remote Backup on MRAWP / Page 26)

2.1.6.3 Auto Report Setup and Remote Backup on MRWP



To create the backup, full access to the system is required. Please inform the customer and request full access. The time required for full access mode is about 30 minutes.

Setup of the CSSI Auto Report at the MRWP



The setup of the CSSI autoreport is only necessary if this function is not available on the system. Please check if the CSSI autoreports are enabled.

ACTION	MRWP
■ Login to the system via cRSP (common Remote Service Platform).	☐
■ Enter the valid service key and click OK .	☐
■ Select Auto Report in the service platform.	☐
■ Under Report , select the entry <code>TEMPLATE</code> .	☐
■ Create the new Auto Report with the name <code>CSSI_Trان_N4</code> Use the following values: - Active: Check this box - Name: <code>CSSI_Trان_N4</code> - Target: <code>SRS_mail</code> - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: <code>C:\MedCom\service\transfer*Numaris4*.ar</code> - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: <code>SRS_ftp</code> - Send when: At once	☐

ACTION

MRWP

Fig. 7: CSSI_Tran_N4


- Click **Save**.

☐

- Confirm the popup with **OK**.

☐

Fig. 8: Save OK

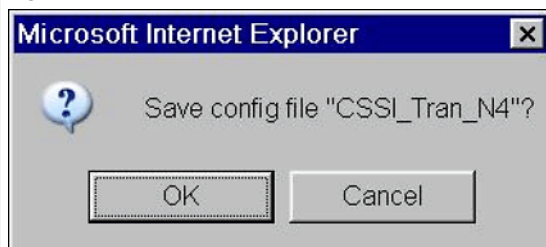

- Continue by clicking **OK** on the popup.

☐

Fig. 9: Success pop-up



- Under **Report**, select the entry **TEMPLATE**.

☐

ACTION	MRWP
■ Create the new Auto Report with the name CSSI_Trان_SecSe Use the following values: - Active: Check this box - Name: CSSI_Trان_SecSe - Target: SRS_mail - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: C:\MedCom\service\transfer*Security-Settings*.ar - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: SRS_ftp - Send when: At once	

Fig. 10: Autoreport CSSI_Trان_SecSe

Siemens Med Service Software - Microsoft Internet Explorer provided by Siemens Medical Solutions

Auto Report

HomeEvent LogHelp

REPORT

TEMPLATE

Active: ☒

Name:

Begin:

12

November

2010

Target:

End:

31

December

2030

Enabling Status:

- ☒ CSSI_Trان_N4
- ☐ IB_Sample
- ☐ IB_Transfer

Sample source?

Source:

Sample what?

Filename:

Search Pattern:

Mode:

Lines:

1

Sample when?

☒ Continuously

Send what?

Mail:

Command:


FTP target:



Send when?

☒ At once ☐ Daily ☐ Weekly ☐ Monthly

SaveDelete

- Click **Save**.

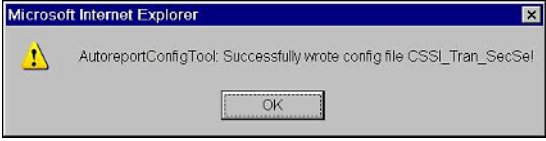ACTION	MRWP
Confirm the popup with OK. Fig. 11: Save OK 	☐
Continue by clicking OK on the popup. Fig. 12: Success pop-up 	☐
Under Report, select the entry TEMPLATE.	☐
- Create the new Auto Report with the name CSSI_Tran_SW02 Use the following values: - Active: Check this box - Name: CSSI_Tran_SW02 - Target: SRS_mail - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: C:\MedCom\service\transfer*SW-Settings*.ar - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: SRS_ftp - Send when: At once	☐

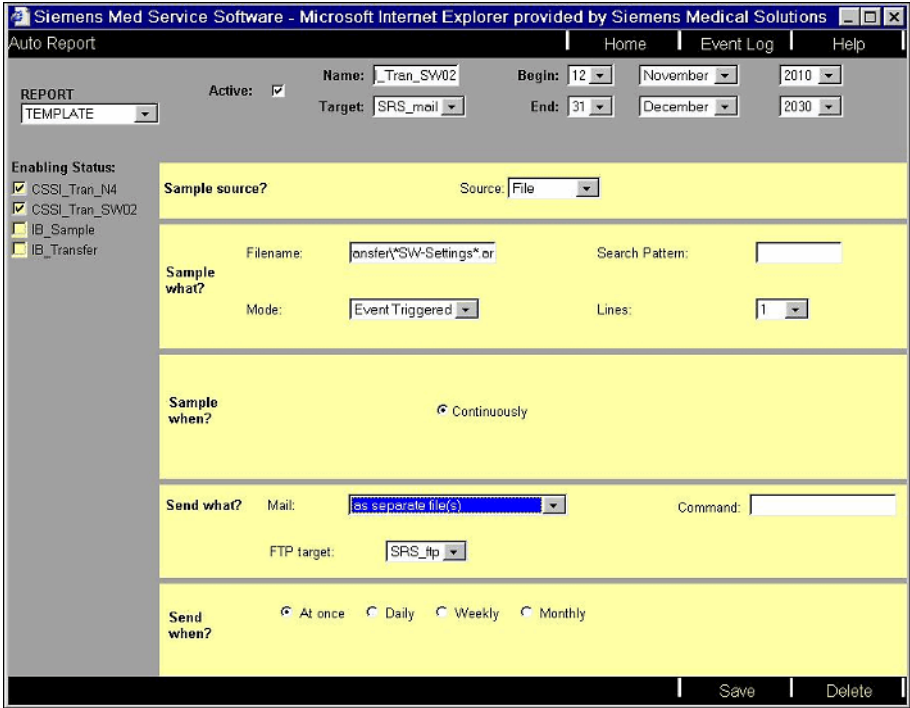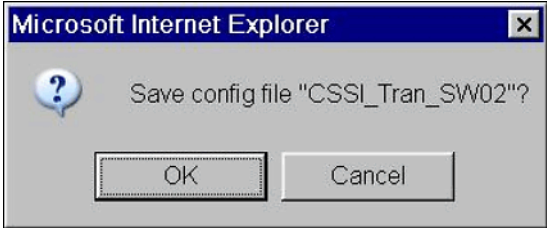ACTION		MRWP
Fig. 13: Autoreport CSSI_Tran_SW02 		
■ Click Save .		☐
■ Confirm the popup with OK .		☐
Fig. 14: Save OK 		
■ Continue by clicking OK on the popup.		☐
Fig. 15: Success pop-up 		
■ Check the newly created Autoreports to see if they are configured correctly and are activated.		☐
■ Click Home to leave the Auto Report section.		☐

Performing the Backup at the MRWP



Before the remote backup can be performed, you must inform the customer about the duration of the backup and the requirement for full access. Also ask the customer to keep the system running in standby mode for the backup transfer.

ACTION	MRWP
Select Backup&Restore in the service platform.	☐
- Perform the backup in local service. The backup has to be performed for the following packages: Numaris4, SW settings02, and Security Settings. Select Backup in the Command menu, select [-E-] FIXED in the Drives menu. - Select SW Settings02 in the Package menu, start the backup by clicking OK. - Select Security Settings in the Package menu, start the backup by clicking OK. - Select Numaris4 in the Package menu, start the backup by clicking OK. Check for errors using the log window.	☐
Click Home to leave Backup&Restore	☐

Move Backups into the transfer directory MRWP

ACTION	MRWP
Select Utilities in the service platform.	☐
Under Source, select the entry Escape to OS in the drop-down menu	☐
Under Command, keep the default entry NT Command Interpreter.	☐
Under Parameters, type in the following: <code>move /y E:*.ar C:\MedCom\service\transfer\</code>	☐

ACTION	MRWP
Click Go. The three backups that have been created are now moved to the directory C:\Med-Com\service\transfer. Fig. 16: Move output	☐
Click Home to leave Utilities	☐
Check the Auto Reports to see whether they are active or not (see (→ Checking active CSSI Auto Reports / Page 34))	☐
Close the service UI in order to close the connection.	☐

2.1.6.4 Auto Report Setup and Remote Backup on MRAWP

Setup of the CSSI Auto Report at the MRAWP

The setup of the CSSI autoreport is only necessary if this function is not available on the system. Please check if the CSSI autoreports are enabled.

ACTION	MRAWP
Login to the system via cRSP (common Remote Service Platform).	☐
Select Auto Report in the service platform.	☐
Under Report, select the entry <code>TEMPLATE</code>.	☐

ACTION	MRAWP
■ Create the new Auto Report with the name CSSI_Trان_N4 Use the following values: - Active: Check this box - Name: CSSI_Trان_N4 - Target: SRS_mail - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: C:\MedCom\service\transfer*Numaris4*.ar - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: SRS_ftp - Send when: At once	☐

Fig. 17: CSSI_Trان_N4

■ Click Save.	☐
--	--------------------------

ACTION	MRAWP
Confirm the popup with OK. Fig. 18: Save OK 	☐
Continue by clicking OK on the popup. Fig. 19: Success pop-up 	☐
Under Report, select the entry <code>TEMPLATE</code>.	☐
- Create the new Auto Report with the name <code>CSSI_Tran_SecSe</code> Use the following values: - Active: Check this box - Name: <code>CSSI_Tran_SecSe</code> - Target: <code>SRS_mail</code> - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: <code>C:\MedCom\service\transfer*Security-Settings*.ar</code> - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: <code>SRS_ftp</code> - Send when: At once	☐

ACTION	MRAWP
---------------	--------------

Fig. 20: Autoreport CSSI_Tran_SecSe


- Click **Save**.

☐

- Confirm the popup with **OK**.

☐


Fig. 21: Save OK

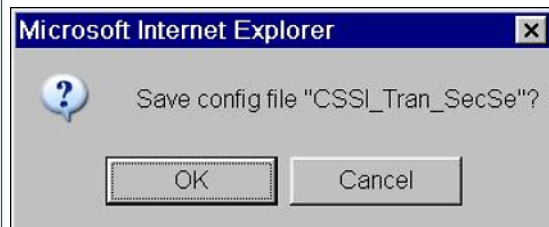

- Continue by clicking **OK** on the popup.

☐

Fig. 22: Success pop-up

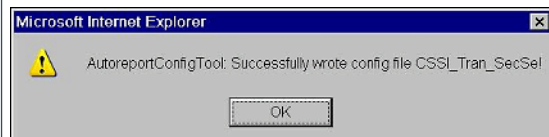


- Under **Report**, select the entry **TEMPLATE**.

☐

ACTION	MRAWP
■ Create the new Auto Report with the name CSSI_Trان_SW02 Use the following values: - Active: Check this box - Name: CSSI_Trان_SW02 - Target: SRS_mail - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: C:\MedCom\service\transfer*SW-Settings*.ar - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: SRS_ftp - Send when: At once	

Fig. 23: Autoreport CSSI_Trان_SW02

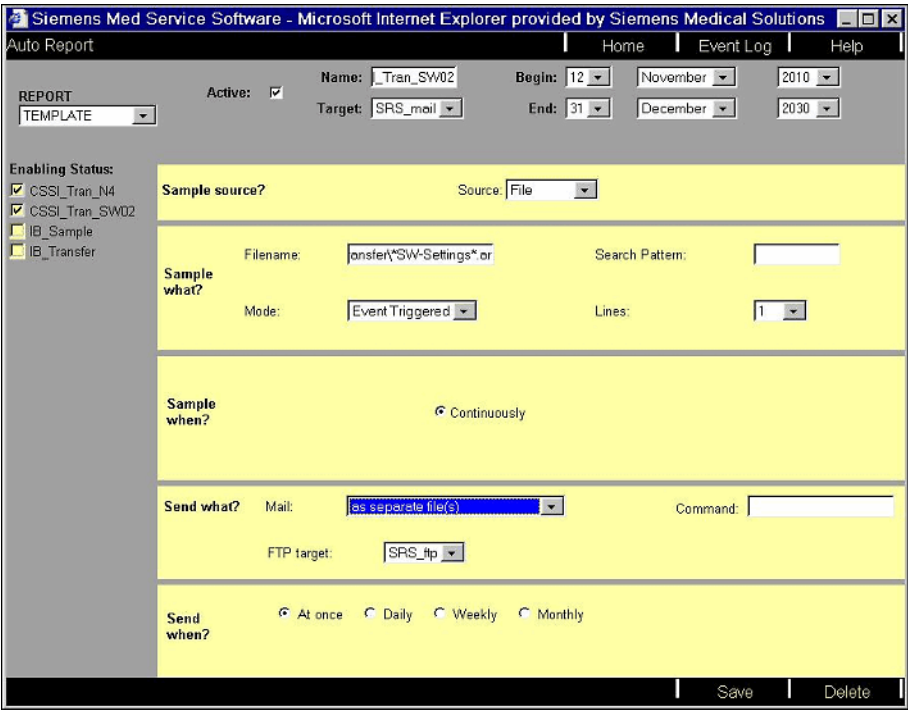

■ Click Save.	
--	------------------------




ACTION	MRAWP
Confirm the popup with OK. <i>Fig. 24: Save OK</i> 	☐
Continue by clicking OK on the popup. <i>Fig. 25: Success pop-up</i> 	☐
Check the newly created Autoreports to see if they are configured correctly and are activated.	☐
Click Home to leave the Auto Report section.	☐

Performing the Backup at the MRAWP



Before the remote backup can be performed, you must inform the customer about the duration of the backup and the requirement for full access. Also ask the customer to keep the system running in standby mode for the backup transfer.

ACTION	MRAWP
Select Reports in the service platform.	☐

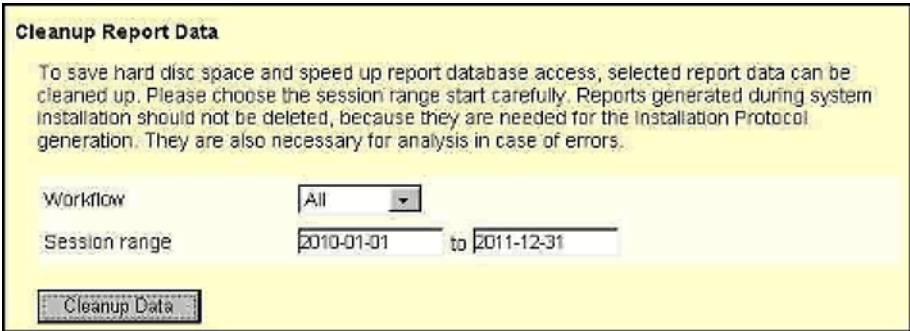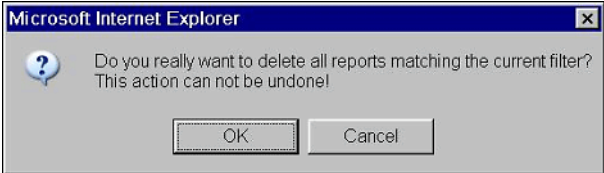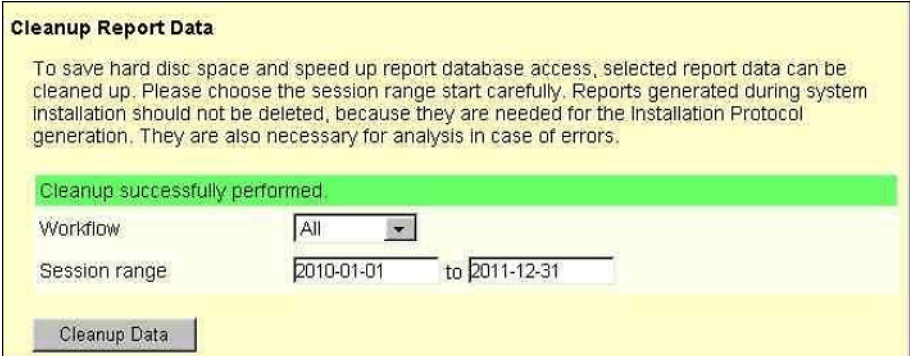ACTION	MRAWP
■ Purge service-related files (test tool files and tune-up, QA files older than one year). If the files are not purged, the size of the master backup created will be huge. - Select Reports -> Clean up - Under Cleanup Report Data select Workflow All - Session range 2010-01-01 to (today - 1 year) e.g. 2011-07-01 - You can also click in the fields and a calendar popup appears. <i>Fig. 26: Data Cleanup</i>  (1) Start date (2) End date - Press the Cleanup Data button. <i>Fig. 27: Confirmation popup</i>  Confirm the popup with OK. - The following shows a successful file purge. <i>Fig. 28: Successful cleanup message</i> 	<input data-bbox="1356 280 1396 324" type="checkbox"/>
■ Click EXIT to close the Report Files popup window	<input data-bbox="1356 1825 1396 1870" type="checkbox"/>
■ Select Backup&Restore in the service platform.	<input data-bbox="1356 1881 1396 1926" type="checkbox"/>

ACTION	MRAWP
- Perform the backup in local service. The backup has to be performed for the following packages: Numaris4, SW settings02, and Security Settings. Select Backup in the Command menu, select [-E-] FIXED in the Drives menu. - Select SW Settings02 in the Package menu, start the backup by clicking OK. - Select Security Settings in the Package menu, start the backup by clicking OK. - Select Numaris4 in the Package menu, start the backup by clicking OK. The MRConvertQuestion popup is not displayed when backup is performed via remote login. The duration of the Numaris4 backup is about 15-20 min. depending on the amount of data. Check for errors using the log window.	☐
Click Home to leave Backup&Restore	☐

Move Backups into the transfer directory

ACTION	MRAWP
Select Utilities in the service platform.	☐
Under Source, select the entry Escape to OS in the drop-down menu	☐
Under Command, keep the default entry NT Command Interpreter.	☐
Under Parameters, type in the following: <pre>move /y E:*.ar C:\MedCom\service\transfer\</pre>	☐

ACTION	MRAWP
■ Click Go. The three backups that have been created are now moved to the directory C:\Med-Com\service\transfer. Fig. 29: Move output	☐
■ Click Home to leave Utilities	☐

2.1.6.5 Checking active CSSI Auto Reports


ACTION	MRAWP/MRWP
■ Select File & Image Tools in the service platform.	☐
■ Under Source, select the entry Transfer List from the drop-down menu.	☐
■ Under Action, select the entry Show/Delete in the drop-down menu.	☐
■ Now the Transfer List is shown. Check to see if the FTP transfer for the backup files has been started. To do this, select the FTP entries in the transfer manager list.	☐

ACTION

MRAWP/MRWP


- For the Numaris4 backup file.


Fig. 30: Transfer List active autoreport Numaris4

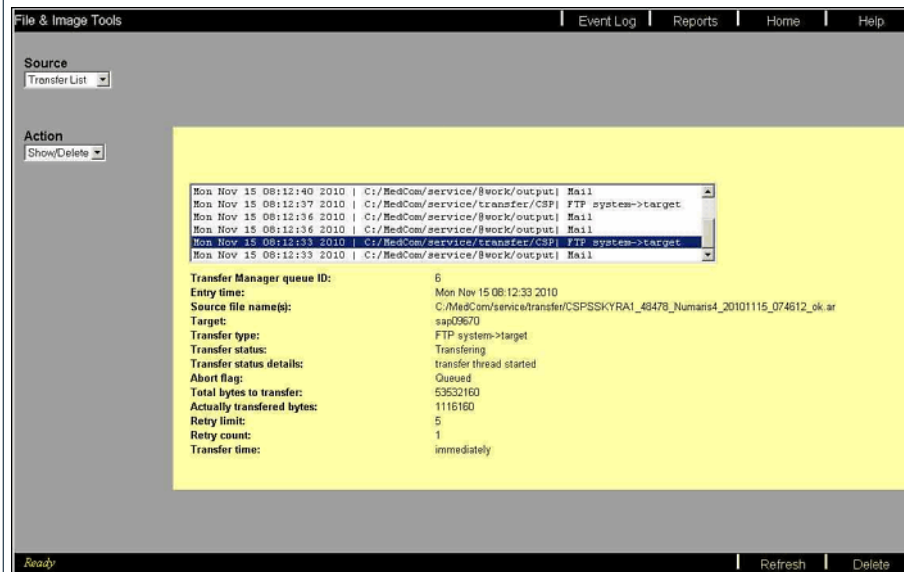

- For the Software02 backup file.

Fig. 31: Transfer List active autoreport Software02

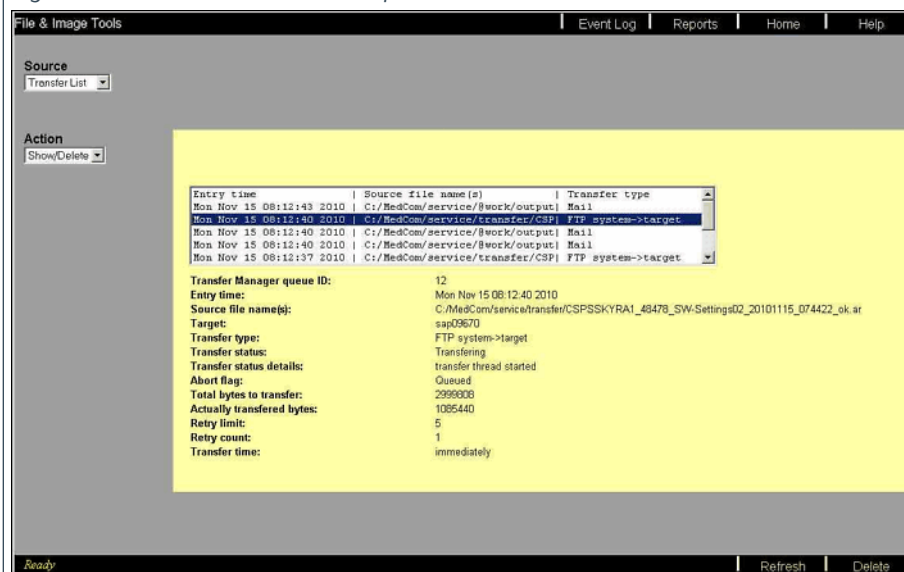

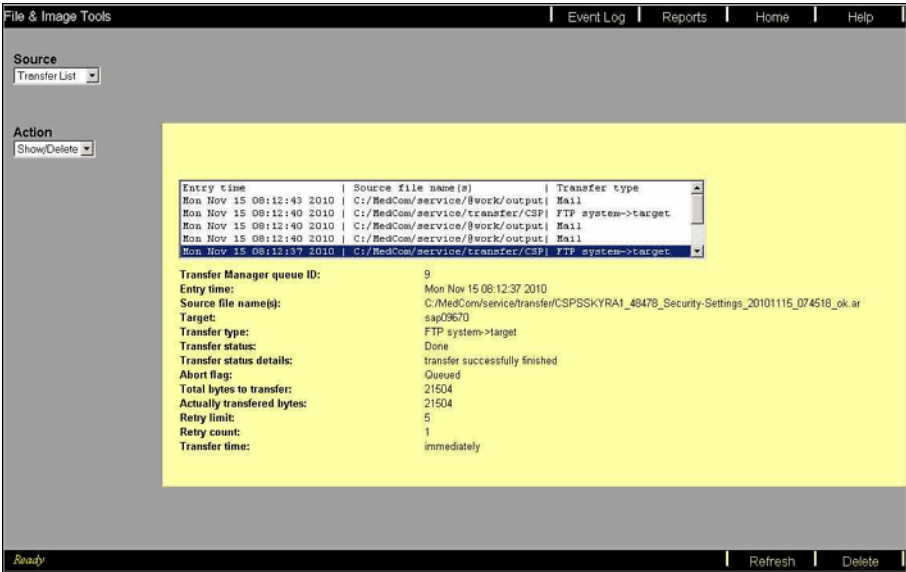ACTION	MRAWP/MRWP
■ For the Security Settings backup file. Fig. 32: Transfer List active autoreport SecuritySettings 	☐
■ Click Home to leave Utilities	☐
■ Close the local service UI.	☐




Please ask the customer not to switch off the overall system. The system should be switched to standby.

3.1 Final Work Steps for EVOLVE from VD13A to VE11D

3.1.1 Overview

- Evolve Installation

Fig. 33: Flowchart - Evolve final steps

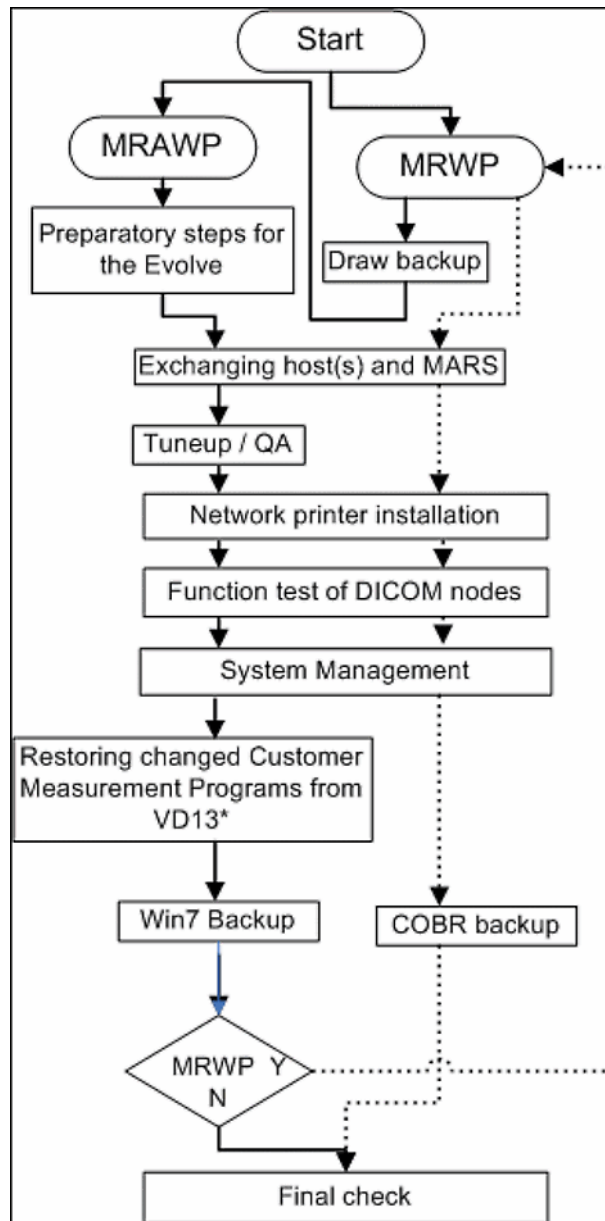

3.1.2 MRWP (MRSC) Preparatory Steps for EVOLVE Update at the MRWP with the Old Software Version



If you have a satellite console with a shared database connected to the main console, perform the following steps first.

If no MRWP is present, skip the following steps and start at (→ MRAWP (MRC) Preparatory Steps for EVOLVE Update at the MRAWP with the Old Software Version / Page 39)



Please check if any local or network printer is configured on MRWP and clarify with the customer representative the availability of installation material (Software, Installation Manual, etc.) if necessary. Please note the configuration data.

3.1.2.1 Performing a Backup on the MRWP



If required (configuration is changed) the step (→ Collecting MRWP Configuration Data from Local Service / Page 12) can be done too.

ACTION	MRWP
Label an empty CD-R medium with Backup VD13A MRWP (actual date)	☐
- Perform the backup in local service. Start local service, select backup under Backup&Restore, insert the labeled CD-R medium into the [CD-R] drive (S:\), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu. - Select Master-Backup in the Package menu, start the backup by clicking OK. Check for errors using the log window.	☐
Remove the labeled CD-R media from the [CD-R] drive (S:\).	☐
Click Home.	☐

3.1.2.2 MNP Deregistration on the MRWP



In order to avoid problems with the MNP registration with the new Numaris version, it is recommended to perform the deregistration with the old Numaris version.

ACTION	MRWP
In the local service UI, select Configuration.	☐
Under System Mgmt. Agent Controls, set the following configuration: Remove the mark from the checkboxes "Go Active" and "Go Online".	☐
Click Go.	☐

ACTION	MRWP
The deregistration of the MNP is done. This will take about 5 minutes.	
■ Click Finish .	☐
■ Click Home and close the local service UI.	☐
■ Shut down the MRWP.	☐
■ Switch off the MRWP.	☐
■ Remove all cables from the MRWP.	☐
■ Continue with (→ MRAWP (MRC) Preparatory Steps for EVOLVE Update at the MRAWP with the Old Software Version / Page 39).	☐

3.1.3 MRAWP (MRC) Preparatory Steps for EVOLVE Update at the MRAWP with the Old Software Version



If you have a satellite console with a shared database connected to the main console, start with the action item:

(→ MRWP (MRSC) Preparatory Steps for EVOLVE Update at the MRWP with the Old Software Version / Page 38)

This has to be done first. If this step has already been completed on the MRWP, continue with the MRAWP.



Please check if any local or network printer is configured on MRAWP and clarify with the customer representative the availability of installation material (Software, Installation Manual, etc.) if necessary. Please note the configuration data.

3.1.3.1 Performing a Backup on the MRAWP



If required (configuration is changed) the step (→ Collecting MRAWP Configuration Data from Local Service / Page 14) can be done too.

ACTION	MRAWP
■ Label an empty CD-R medium with Backup VD13A MRAWP (actual date)	☐
■ Perform the backup in local service. Start local service, select backup under Backup&Restore , insert the labeled CD-R medium into the [CD-R] drive (S:I), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu. - Select Master-Backup in the Package menu, start the backup by clicking OK . Check for errors using the log window.	☐

ACTION	MRAWP
Remove the labeled CD-R media from the [CD-R] drive (S:\).	☐
Click Home.	☐

3.1.3.2 MNP Deregistration on the MRAWP



In order to avoid problems with the MNP registration with the new Numaris version, it is recommended to perform the deregistration with the old Numaris version.

ACTION	MRAWP
In the local service UI, select Configuration.	☐
Under System Mgmt. Agent Controls, set the following configuration: Remove the mark from the checkboxes "Go Active" and "Go Online".	☐
Click Go.	☐
The deregistration of the MNP is done. This will take about 5 minutes.	
Click Finish.	☐
Click Home and close the local service UI.	☐
Close the local service UI.	☐

3.1.3.3 Manual backup of voice commands on MRAWP



Due to a bug in the backup program, it is recommended to copy the voice command files to an external medium and then to restore these in the new software version.

Please check the external medium first with the virus scanner of the service PC. Ensure that the external medium used is free of viruses.

ACTION	MRAWP
- Enable the advanced user - Run the Advanced User program from the Start menu (press <CTRL> + <ESC> to access the Start menu) - From the Start menu, select All Programs -> Advanced User. - In the popup displayed, enter the advanced user password and press the OK button - Confirm Advanced User enabled popup with OK.	☐
Connect a USB medium.	☐

ACTION	MRAWP
■ Start Windows Explorer: - Press <CTRL> + <ESC> to access the Start menu - In the Start menu select All Programs ⇒ Accessories ⇒ Windows Explorer Windows Explorer starts.	☐
■ Navigate to the directory C:\MedCom\MriSiteData	☐
■ Copy the subdirectory PIA and its contents to the USB medium. The subdirectory PIA contains the voice command files.	☐
■ Close Windows Explorer.	☐
■ Remove the USB medium.	☐
■ Perform a shutdown of the system.	☐
■ Switch off the system on the alarm box.	☐
■ Remove all cables from the MRAWP.	☐

3.1.4 Replacing the MRAWP/MRWP and MARS



Adhere to the C-MOS guidelines when handling components.

3.1.4.1 Required Tools and Auxiliary Materials

- Anti-static wristband and grounded tabletop pad
(TDK equipment, Type 8005 from 3M Deutschland GmbH)

3.1.4.2 Replacing the MRWP

The following procedure has to be performed:

- Reuse the old USB dongle for the update
 - Remove the USB dongle from the old MRAWP.
 - Plug the old USB dongle from the old into the new HPz440.

- Reconnect all cables on the new MRWP (HPz440)
All other connections see sticker on "Labelling of external Connectors"

Fig. 34: Network connection of the MRWP



(1) Hospital Network


3.1.4.3 Replacing the MRAWP

The following procedure has to be performed:

- Reuse the old USB dongle for the update.
 - Remove the USB dongle from the old MRAWP.
 - Plug the old USB dongle from the old into the new HPz440.

- Reconnect all cables on the new MRAWP (HPz440).

Fig. 35: Network connection of the MRAWP

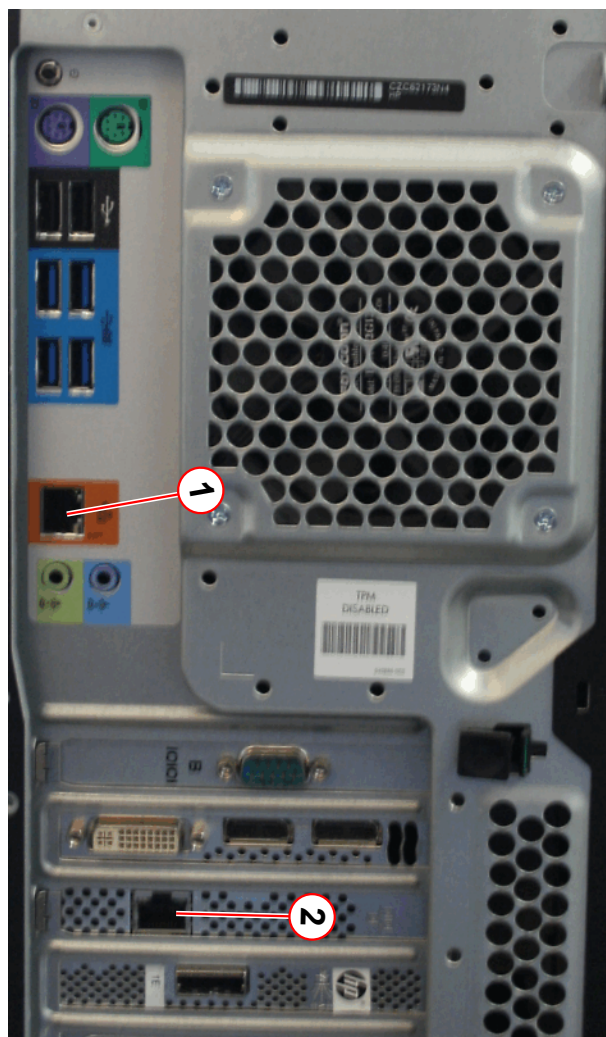

- (1) Hospital Network (external)
- (2) Measurement Network (internal)

All other connections see sticker on "Labeling of external Connectors"

Fig. 36: MOXA Box




MRAWP with new 2-port MOXA USB device:

- Remove the old 4-port connector/distribution box of the old MOXA card.
- Connect the USB cable of the 2-port MOXA USB device according the lable on the HPz440 host.
- Assign and connect the device(s) to the new 2-port MOXA card according to the following table. The port number is assigned on the connectors of the 4x-port adapter/connector cable:

Old cable connector port	Assigned device	Assigned COM port	Port number on the MOXA distribution box
P4	Alarm Box	COM9	Port 2
P3	ACC (MPS)	COM6	Port 1

- Use the delivered adapter cable 9 pin to 25 pin to connect the 25 pin cable from the alarmbox to the new USB Moxa device.
Use a cable tie to ensure a stable connection.
- If the extension cable is needed, remove the attachment screws of the connector (pull out and turn counter-clockwise) (→ 1/ Fig. 37 Page 45) at the female connector of the extension cable. Subsequently, secure the attached hex-headed extension bolts on the front side of the connector(→ 2/ Fig. 37 Page 45).


Fig. 37: Female connector of the extension cable

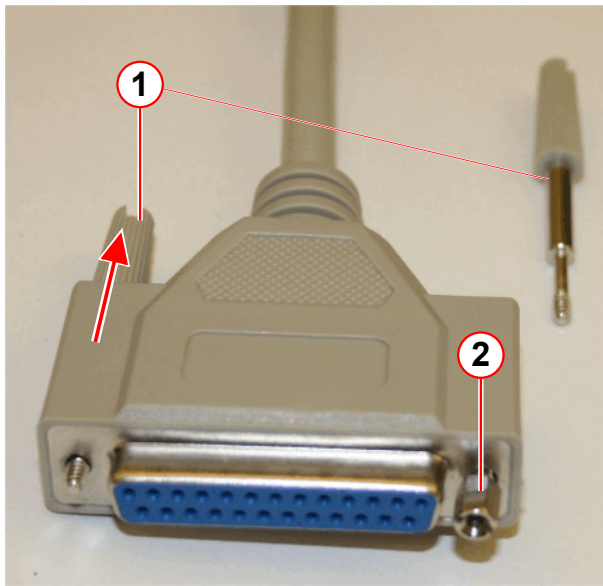

- (1) Attachment screw for connector
- (2) Hex-headed extension bolt

- (→ Fig. 38 Page 45) shows when the extension cable is connected with the 4x-adapt-er/connector cable of the new 4-port MOXA card.

Fig. 38: Connected extension cable

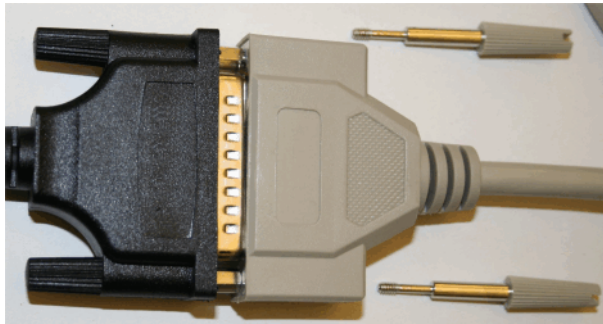

- Connect the external device with the new USB MOXA device. If needed, replace the screws on the other connector of the extension cable the same way to ensure a stable connection.

3.1.4.4 Replacing the MARS



The new delivered MaRS0S.7 18/32RX part.#11267844 supports both RX channel configurations (18 and 32 Channels).

The following procedure has to be performed:

- Remove all cables on the old installed MARS.
- Remove the old MARS from the system cabinet.
- Insert the new MARS supplied into the system cabinet.

- Reconnect all cables on the MARS.



Be careful to connect the fiber optic cables to the OPEC cards. Check the correct PCI-e slot location of the OPEC cards. To do this, a label can be found on the new MARS.

3.1.5 MRAWP (MRC) System Start-up



With E11D the new Password Handling in syngo is introduced. Further information see (→ New Password Handling in syngo / Page 94)



In some cases the following popups might be displayed. They can be ignored.

Fig. 39: Reboot popup


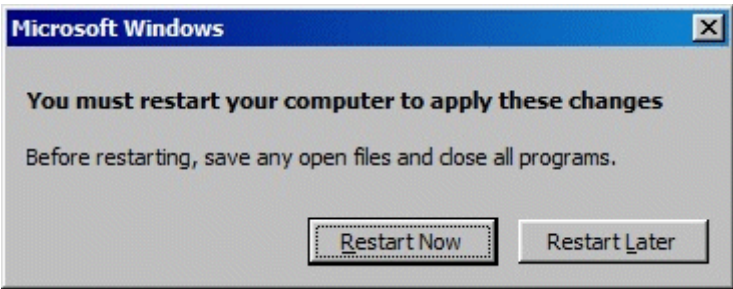


Fig. 40: Network location popup

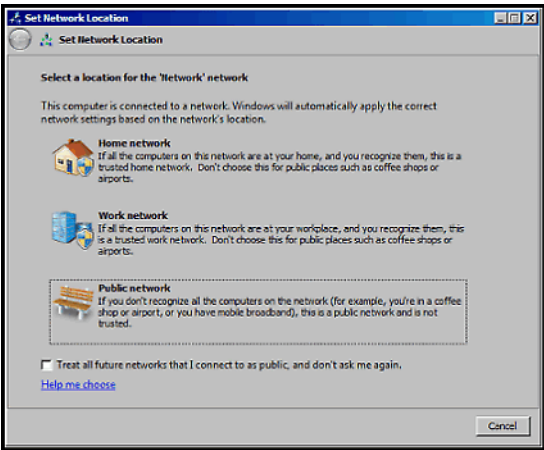


ACTION	MRAWP
■ Switch on the MR system at the alarm box.	☐
■ Wait until the main UI is displayed and wait until the measurement system is ready. The firmware update time duration is about 30 minutes.	☐


3.1.6 MRAWP (MRC) New secure Service Key Upload

- 

With VE11D a new secure Service Key is introduced. More information can be found here (→ Service key handling with new Secure Service Key / Page 89)
- 

The service key with copy and from the “*ServiceKey.dat” file. The file can be found in the following subdirectory C:\temp\license. Please delete this directory permanent after the update and store the service key file on a different location (e.g on the USB hardisk which is used for the COBR image).

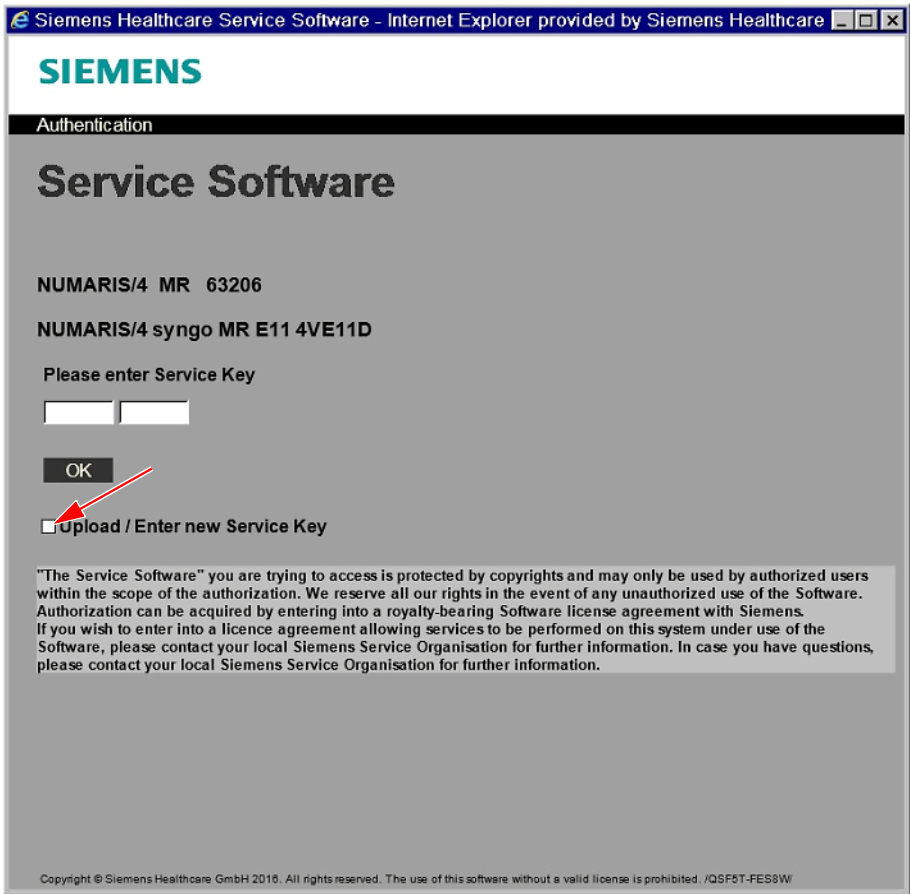ACTION	MRAWP
■ Select Options ⇒ Service ⇒ Local Service	☐
■ Select the checkbox “Upload/Enter new Service Key” to enable the upload <i>Fig. 41: Mark checkbox Upload</i> 	☐

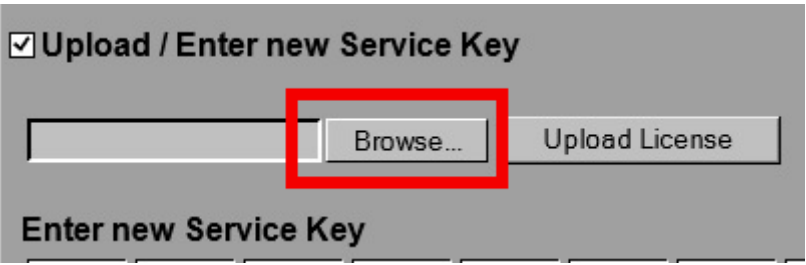ACTION	MRAWP
Select the button Browse and navigate to the location of the <code>*ServiceKey.dat</code> file and select the file. <i>Fig. 42: Browse button</i>  The screenshot shows a dialog box titled 'Upload / Enter new Service Key'. It contains a text input field, a 'Browse...' button (highlighted with a red rectangle), and an 'Upload License' button. Below the input field is a section titled 'Enter new Service Key'.	☐
Select the button Upload License.	☐
Now the service key is uploaded and a login to the service UI is done.	☐


3.1.7 MRAWP (MRC) Configuration checks





By default the advanced user is disabled. For the configuration check the advanced user must be enabled (see (→ MRAWP (MRC) Setup Advanced User Login Function / Page 65).

ACTION	MRAWP
Select Options ⇒ Service ⇒ Local Service.	☐
Enter the valid service key and click OK.	☐
Select Configuration in the service platform.	☐

ACTION	MRAWP
■ Select Monitor Type under Local Host and check the settings as follows: - Monitor Type: Color and Siemens (True Color) - Multi-Monitor settings Only "Start syngo in single monitor" is marked Fig. 43: Monitor settings SESO 	☐
■ Click Save if changes have been made and Finish	☐

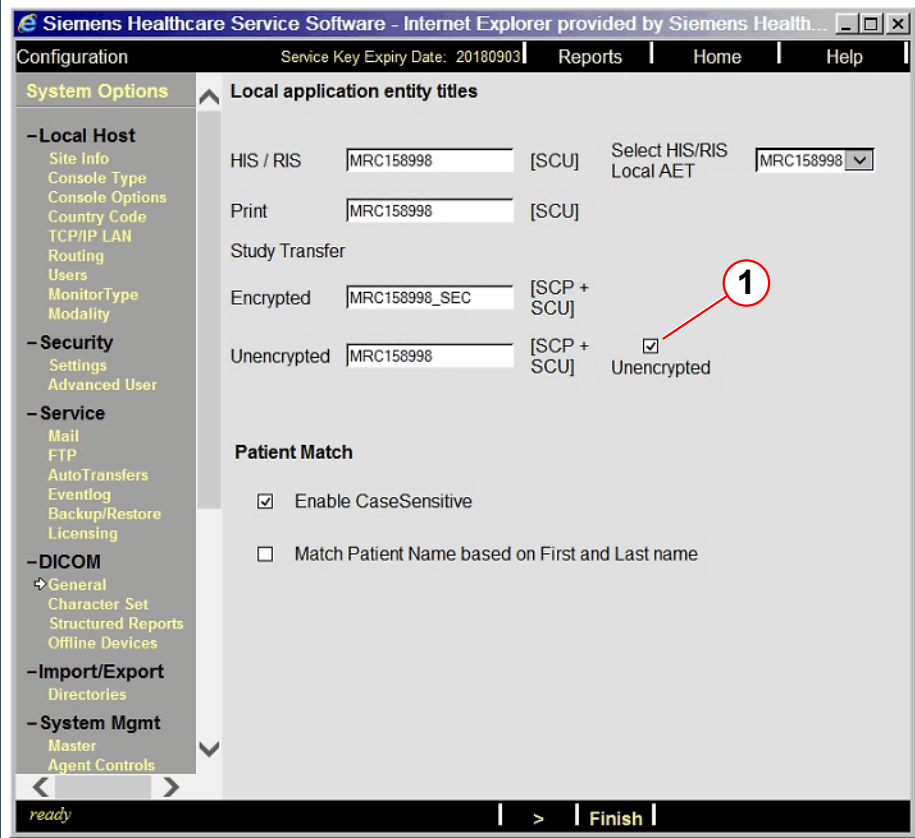ACTION		MRAWP
■ Under DICOM, select General and check if the checkbox Unencrypted is marked. Fig. 44: Example from SIM Console  (1) Unencrypted checkbox must be marked		☐
■ Click Save if changes have been made and Finish		☐
■ Under DICOM, select Network Nodes and press under Host Properties >.		☐
■ On the next page, under General Node Properties in the drop-down menu Logical Name, select each available entry. For each entry check the Connection Mode. This must be set to Unencrypted. Fig. 45: DICOM connection mode  Click Save if changes have been made for each entry. All configured DICOM nodes from the former software version can be found under C:\MedCom\MriSiteData\report.pdf or C:\MedCom\MriSiteData\BaRe_HW_Info\MRC\report.html If there are some missing, please carry out the configuration again.		☐
■ If all entries have been checked, click Finish.		☐
■ Under DICOM, select Print Devices and press under Host Properties >.		☐

ACTION	MRAWP
On the next page, under General Node Properties in the drop-down menu Logical Name, select each available entry. For each entry check the Connection Mode. This must be set to Unencrypted. <i>Fig. 46: DICOM connection mode</i>  Click Save if changes have been made for each entry.	☐
If all entries have been checked, click Finish.	☐
Under DICOM, select HIS/RIS Nodes and press under Host Properties >.	☐
On the next page, under General Node Properties in the drop-down menu Logical Name, select each available entry. For each entry check the Connection Mode. This must be set to Unencrypted. <i>Fig. 47: DICOM connection mode</i>  Click Save if changes have been made for each entry.	☐
If all entries have been checked, click Finish.	☐

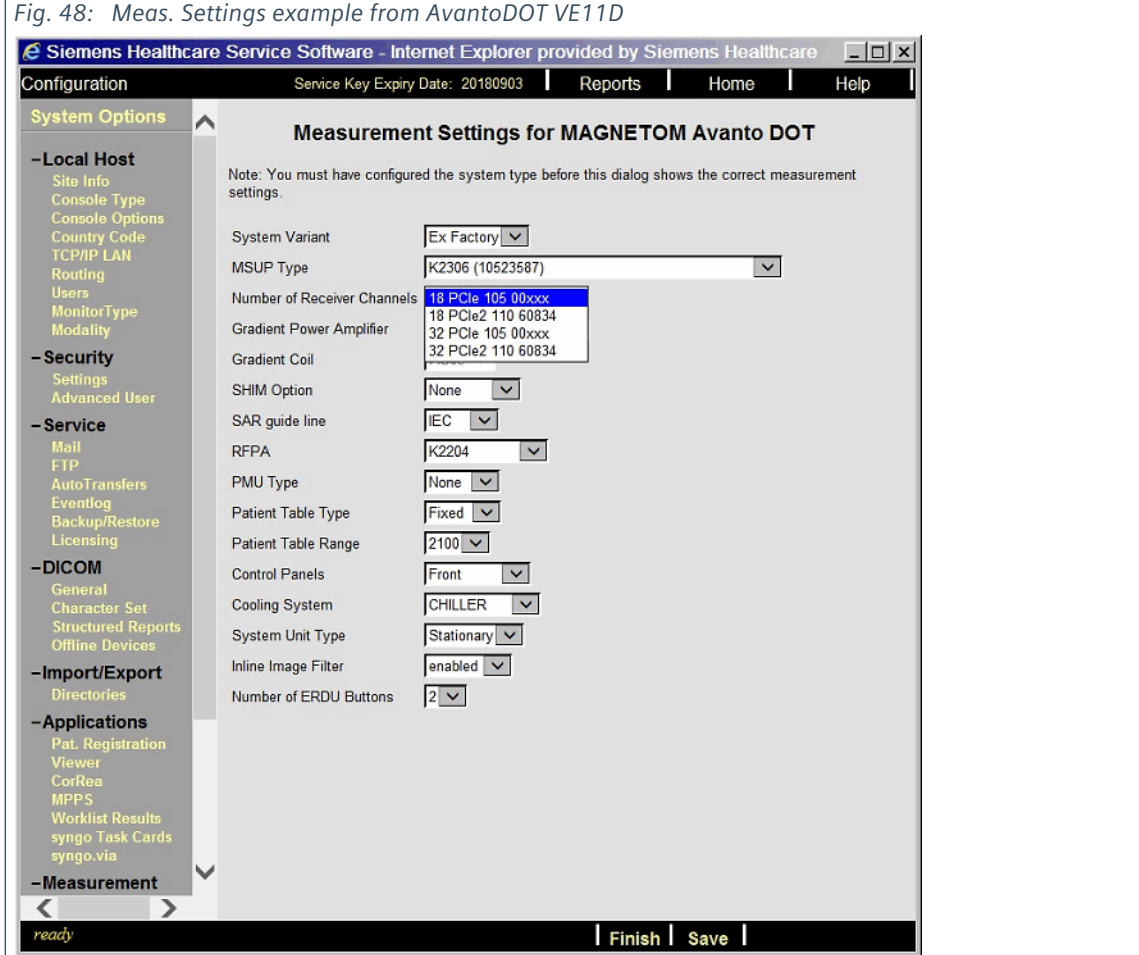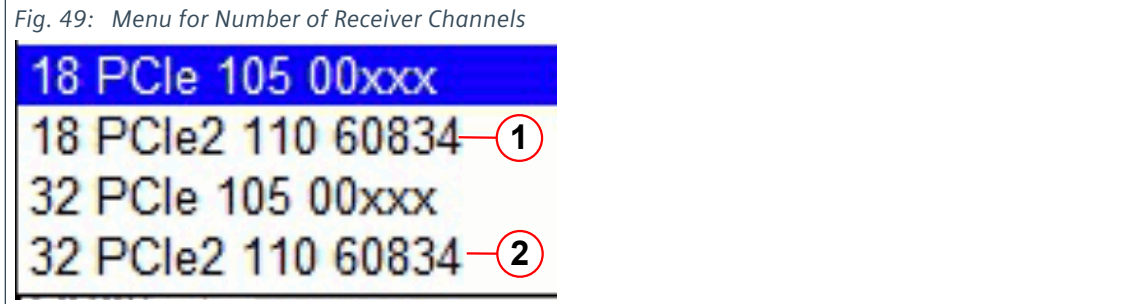ACTION	MRAWP
■ Under Measurement Meas. Settings ,check the settings according to the MR Report from VD13A. C:\MedCom\MriSiteData\BaRe_HW_Info\MRC\MrReport.htm C : \MedCom\MriSiteData\BaRe_HW_Info\MRC\MrReport.htm Fig. 48: Meas. Settings example from AvantoDOT VE11D  With evolve a new MaRS (MaRS0S.7 18/32RX part.# 11267844i) s delivered with new PCIe2-RX16 cards. Check if the correct PCIe2-RX16 entry according the receive channels is sected under Number of Receiver Channels. Fig. 49: Menu for Number of Receiver Channels  (1) Section for new MaRS with 18 Rec. Channels (2) Section for new MaRS with 32 Rec. Channels Click Save if changes have been made for each entry.	☐ ☐
■ Click Finish to continue.	☐


ACTION	MRAWP
Click Home and confirm the requested restart of application with OK.	☐
Wait until the restart is performed and the syngo main UI is up and running.	☐



If a new Network Node, Print device or HIS/RIS Node is configured, be aware that the default setting for the Connection Mode is Unencrypted. Please keep it Unencrypted.

3.1.8 MRAWP (MRC) Tune-up and QA Steps

ACTION	MRAWP
Select Options --> Service --> Local Service.	☐
Enter the last 10 characters from the service key into two shown groups.. <i>Fig. 50: Service Key groups</i>  This can be done either by typing in each character or via copy/paste from the service key file. Press OK to login.	☐

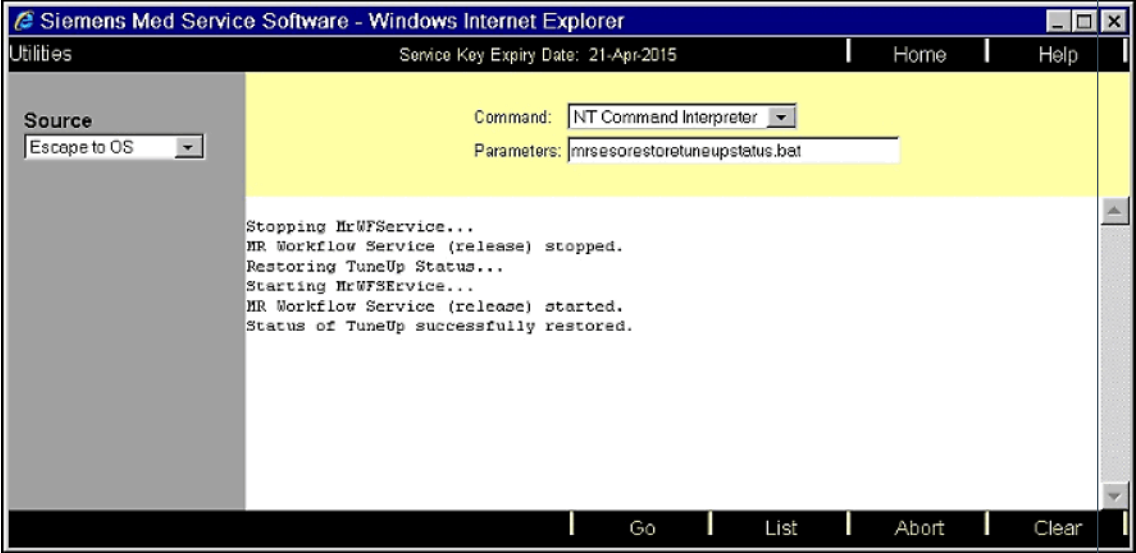ACTION	MRAWP
- Restore the Tuneup status: - Select Utilities ⇒ Escape to OS and type in the command <code>MrSesoRestoreTuneUpStatus.bat</code> and press <ENTER> Fig. 51: SESO output  This step restores the last TuneUp&QA steps. The QA steps for the surface coils are restored too.	☐
Click Home and select Tune Up.	☐
Perform all tune-up steps that have the status "ToDo"	☐
Note: If you have a satellite console MRWP (MRSC), you can now start with the (→ MRWP (MRSC) Startup / Page 67) step in parallel with the tune-up and QA.	
- Perform all QA steps that have the "ToDo" status If all QA steps are set to "Done" perform the following steps: - BC Power Losses Check - BC Image Check - Image Orientation Check	☐



Before closing and deregister the service patient by the Service Software wait until a running measurement (e.g. Tune Up or QA) is stopped (normal end or manual user stop). Otherwise an internal error occurs and the exam task card crashes. If this happens a re-boot is required.

3.1.9 MRAWP (MRC) Surface Coils QA Steps



For all coils used before the software update, it is not necessary to perform a Coil QA. The QA status is restored.

3.1.10 MRAWP (MRC) Network Printer Installation

ACTION	MRAWP
Install the network paper printer according to the software installation manual (this step is not part of the dynamic data) if previously installed on the MRAWP or MRWP.	☐
If applicable, set the installed printer as the default printer. Select Start -> Settings -> Printers and Faxes Select the printer and click the right mouse button. In the menu displayed, select the "Default printer" entry.	☐

3.1.11 MRAWP (MRC) Function Test of the DICOM Nodes

ACTION	MRAWP
Send an image from the Patient Browser to all connected DICOM nodes. The images must be sent successfully	☐
Send a film sheet to all connected DICOM cameras. The film sheet must be sent successfully to the DICOM camera.	☐
Retrieve a test patient from the RIS/HIS system, if available. The patient should be available on the system.	☐



A DICOM connection can now be made secure (encrypted). This is a new configuration parameter. See (→ New feature for DICOM node configuration / Page 100)



The new DICOM Node Configuration Option "Supports Presentation State" must not be enabled for any DICOM storage node. The NUMARIS/4 VE11 does not currently support the "Presentation State" for created MR images.


3.1.12 MRAWP (MRC) System Management





Due to the CSSI process, system management cannot be activated at the factory. System management has to be activated at the customer site.


3.1.12.1 Installation and Activation Steps


ACTION	MRAWP
Select Options --> Service -> Local Service.	☐
Enter the valid service key.	☐
In the local service UI, select Configuration	☐
Check that the settings for service region (Europe, USA, or Asia) under Local Host Site Info are correct. Correct the settings, if necessary, and save the changes. This entry assigns the IP addresses for the RDIAG server and the System Management server.	☐
Under the "Local Host" Routing, type in the IP address of the router or gateway used for SRS. In some cases, this entry is the default gateway that is used for the TCP/IP LAN configuration. See log file <code>C:\MedCom\MriSiteData\ExtBaRe\routing.txt</code> for the routing entry. For the installation of System Management, the correct entry under routing is essential. If this routing entry is not set (empty), you will get an error popup when entering the System Management configuration.	☐
Select System Options. Fig. 52: Configuration 	☐

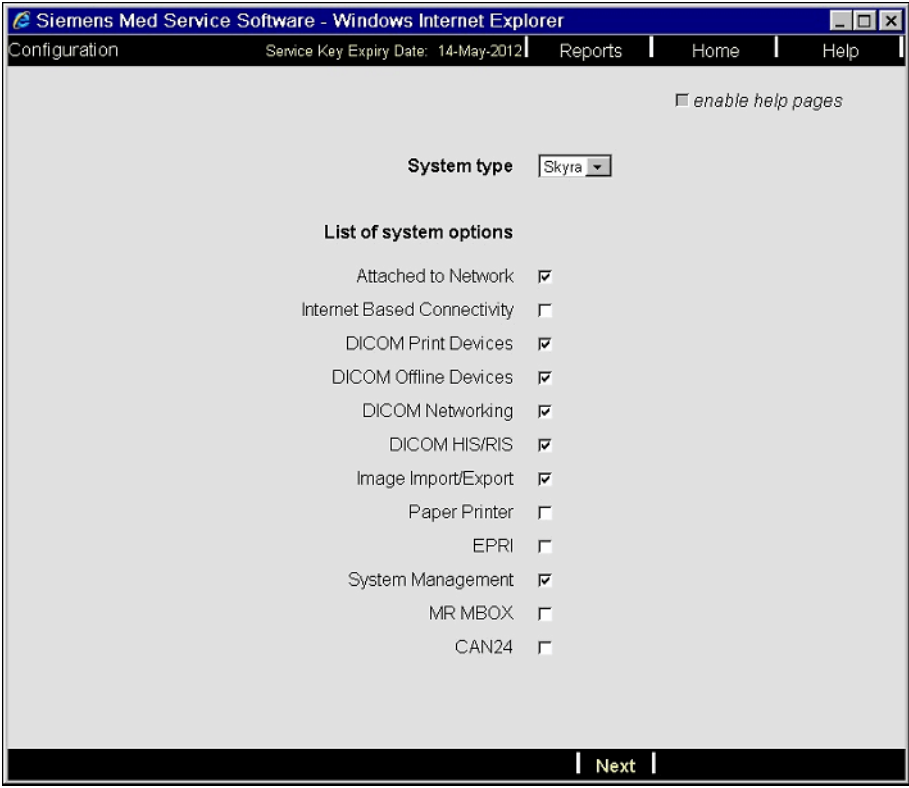ACTION	MRAWP
Under "List of System Options", click the System Management checkbox Fig. 53: System options  The screenshot shows the 'Siemens Med Service Software - Windows Internet Explorer' window. The address bar shows 'Configuration'. The page has a navigation bar with 'Service Key Expiry Date: 14-May-2012', 'Reports', 'Home', and 'Help'. A checkbox labeled 'enable help pages' is in the top right. Below this is a 'System type' dropdown menu currently set to 'Skyra'. The main section is titled 'List of system options' and contains a list of features with checkboxes: 'Attached to Network' (checked), 'Internet Based Connectivity' (unchecked), 'DICOM Print Devices' (checked), 'DICOM Offline Devices' (checked), 'DICOM Networking' (checked), 'DICOM HIS/RIS' (checked), 'Image Import/Export' (checked), 'Paper Printer' (unchecked), 'EPRI' (unchecked), 'System Management' (checked), 'MR MBOX' (unchecked), and 'CAN24' (unchecked). At the bottom of the page is a 'Next' button.	☐
Select Next.	☐

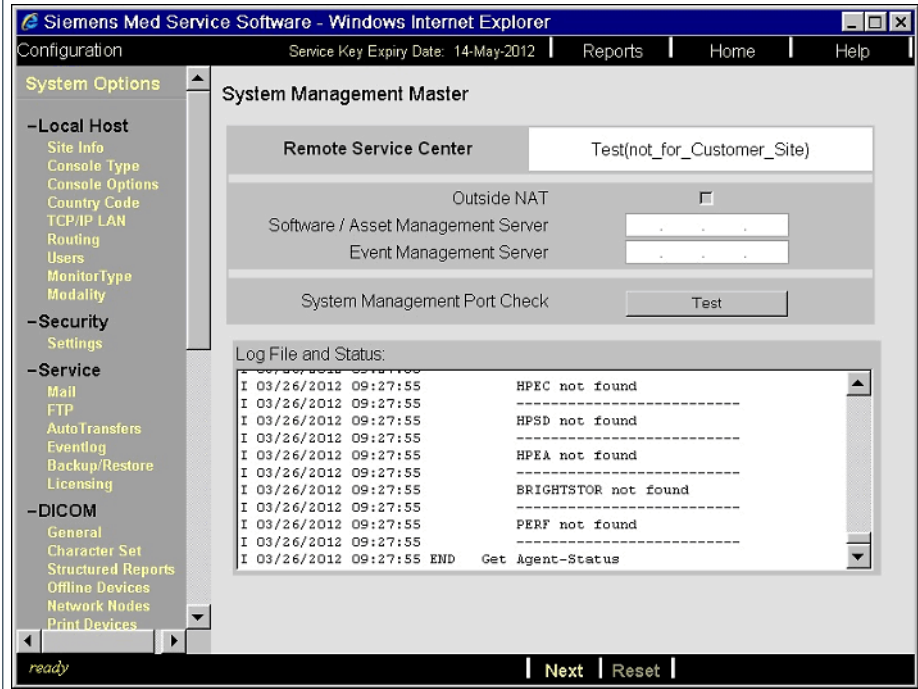ACTION	MRAWP
■ On the System Management Master configuration page Fig. 54: System Management master page  See the report file C:\MedCom\MriSiteData\BaRe_HW_Info\MRC\report.html for the Outside NAT configuration.	☐

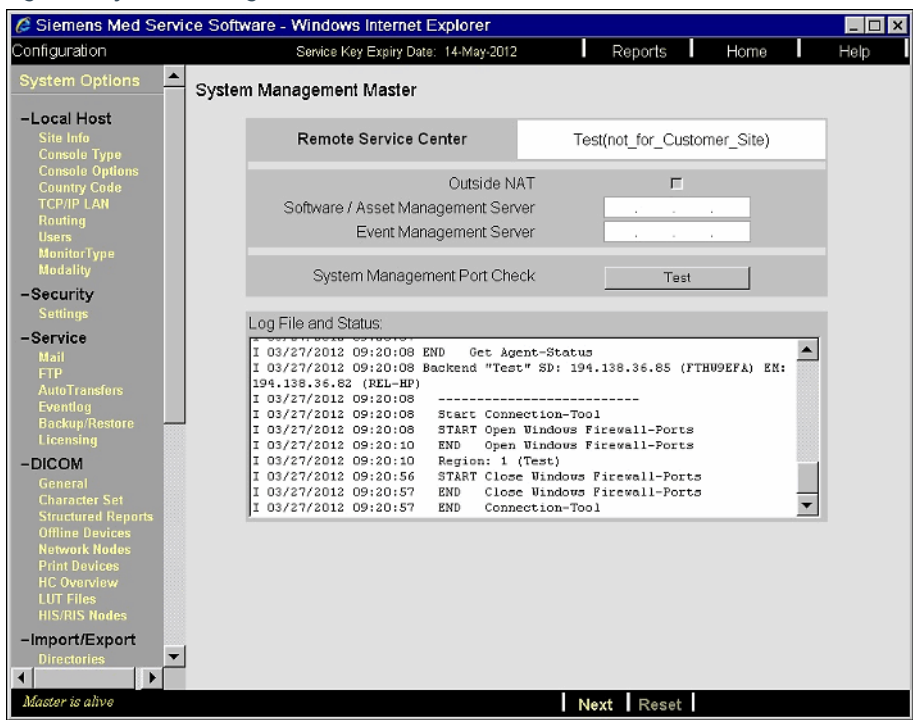ACTION	MRAWP
■ Use the Test button to check whether the System Management server can be accessed. Fig. 55: System Management master can be reached  The log file for the test can be found under the directory: C:\MedCom\service\mwtools\CCT\Logs	☐
In the event of problems with the connection test, contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/ Technical Support Group.	
■ Select Next.	☐

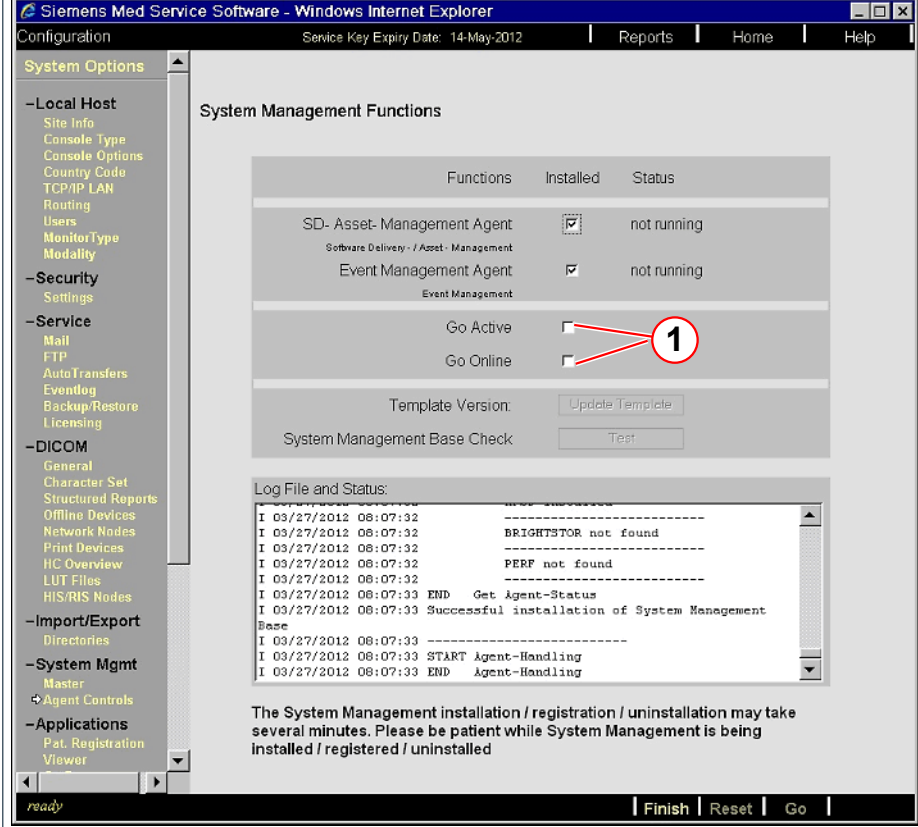ACTION	MRAWP
Under System Mgmt. Agent Controls, set the following configuration: Mark the checkboxes "Go Active" and "Go Online". (The other checkboxes have already been marked). Fig. 56: System Management activation  (1) Mark both checkboxes In the Log File and Status area, the current status of System Management activation/ installation is displayed.	☐ ☐
Click Go. Note: The installation/registration will take approximately 10 to 15 minutes.	☐

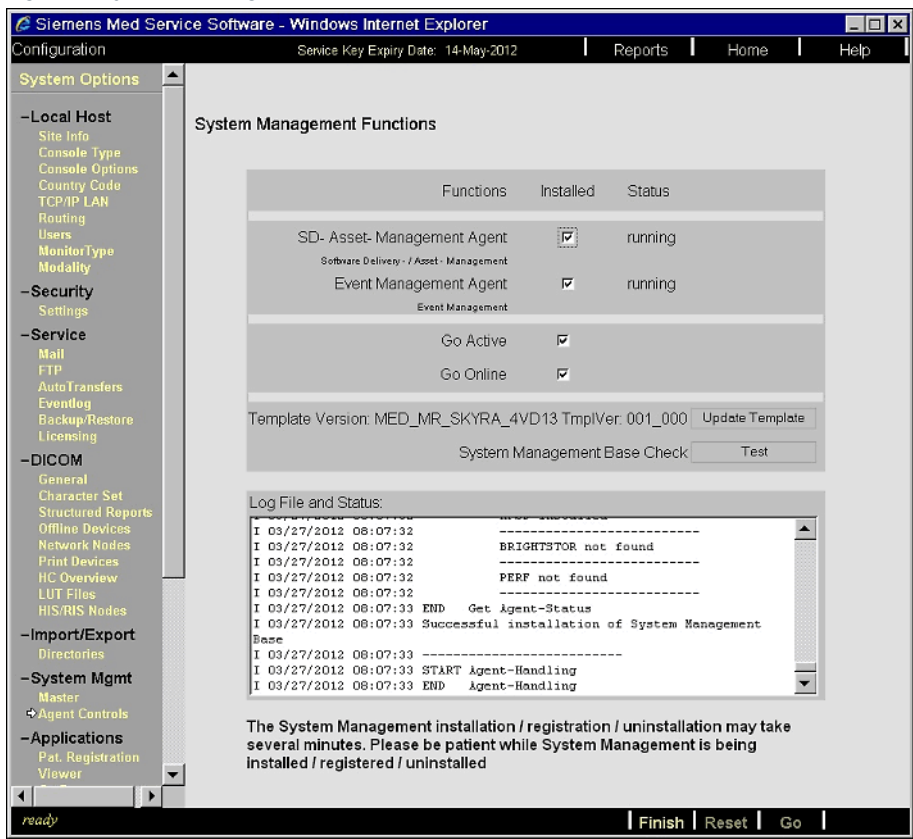ACTION	MRAWP
■ If activation is successful, the following is displayed: <i>Fig. 57: System Management activation OK</i>  Select Finish to continue	☐
■ Click Finish and close the local service UI. In the event of problems (installation/activation fails a second time), contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/Technical Support Group.	☐





For problems installing the Event Management Agent, see (→ System management installation workaround / Page 139)

3.1.13 MRAWP (MRC) Restoring Changed Customer Measurement Programs from VD13A



The following step restores all changed customer protocols. Only these protocols will be converted after the reboot.

ACTION	MRAWP
Insert the MasterBackup CD labeled Master Backup VD13A MRAWP (actual date) created; see (→ Performing a Backup on the MRAWP / Page 39) into the CD/DVD-ROM drive R:	☐
Select Options --> Service --> Local Service.	☐
Enter the valid service key and click OK.	☐
Select Backup&Restore in the service platform.	☐
Select under Command Restore	☐
Under Drives, select [-R-] IF the drive letter not correct, please perform (→ Drive Letter Assignment for the CD/DVD-ROM Drive and/or for the CD/DVD-R Drive / Page 101).	☐
Under Archive, select the MasterBackup and under Groups the entry N4Customer.	☐
Press Go to start the restore. Wait until the restore has finished.	☐
Remove the Master Backup CD labeled Master Backup VD13A MRAWP (actual date) from the CD/DVD-ROM drive R:	☐
Click Home and confirm the requested reboot message popup with OK.	☐

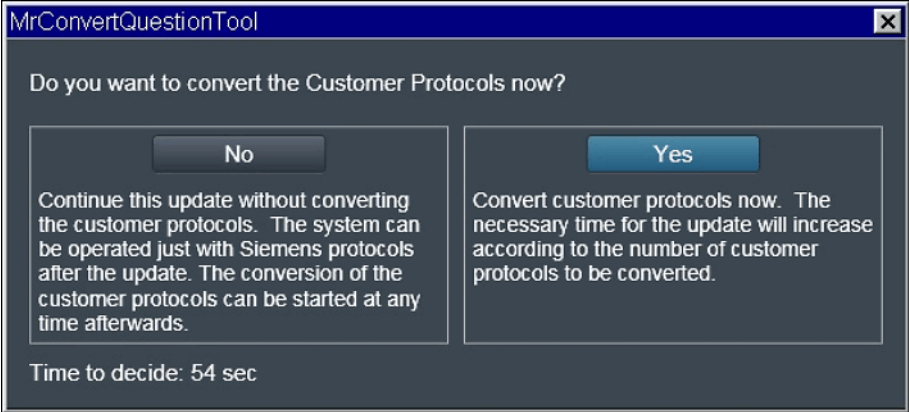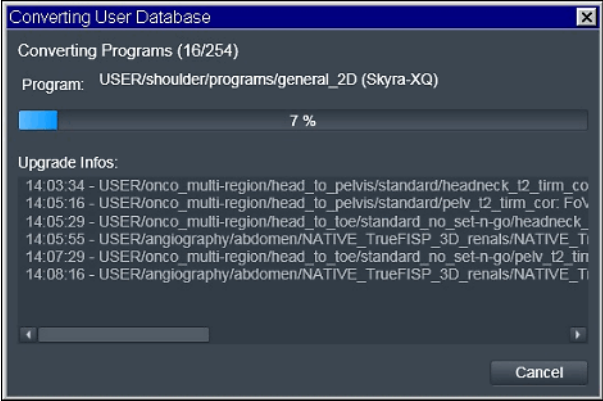ACTION	MRAWP
A reboot is now performed. After the reboot the changed protocols will be converted. Confirm the MRConversionQuestion popup with Yes. Fig. 58: MRConversionQuestion popup  Now the conversion will start. Fig. 59: Convert popup  Please wait until this is finished.	☐
Wait until the main UI is up and the system is ready for measurement.	☐

3.1.14 MRAWP (MRC) Manual Restore of Voice Commands



Due to a bug in the backup program, it is recommended to restore the voice command files from an external medium in the new software version.

Please check the external medium first with the virus scanner of the service PC. Ensure that the external medium used is free of viruses.

ACTION	MRAWP
■ Enable the advanced user - Run the Advanced User program from the Start menu (press <CTRL> + <ESC> to access the Start menu) - From the Start menu, select All Programs -> Advanced User. - In the popup displayed, enter the advanced user password and press the OK button - Confirm Advanced User enabled popup with OK.	☐
■ Connect the USB medium containing the backup of the voice command files.	☐
■ Start Windows Explorer: - Press <CTRL> + <ESC> to access the Start menu - In the Start menu select All Programs ⇒ Accessories ⇒ Windows Explorer Windows Explorer starts.	☐
■ Navigate to the USB medium	☐
■ Copy the directory PIA and its contents from the USB medium to the directory C:\MedCom\MriSiteData	☐
■ Confirm folder replacement popup with Yes	☐
■ In the Copy File popup mark the checkbox "Do this for the next ... conflicts" and select Copy and Replace. The voice command files are now copied into the directory C:\MedCom\MriSiteData\PIA	☐
■ Close Windows Explorer.	☐
■ Remove the USB medium.	☐

3.1.15 MRAWP (MRC) Standby Function Check

Check if the "Standby" status can be activated via the System Manager and restart the MR system as follows:

- Select **System --> Control**, open the tap MR Scanner, and click the Standby button.
- Check to see if the scanner is turned off.

This will take a few minutes due to the automatic shutdown of the imager. If the scanner is not turned off, please check the com port settings, refer to (→ Checking and setting the MOXA - Ports for HPz440 / Page 148)
- Then press the System On button After a while, the information 'Measurement system ready' will appear.

3.1.16 MRAWP (MRC) Setup Advanced User Login Function



By default the Advanced User Login function is disabled. If it needs to be changed, see (→ Advanced User handling / Page 98)



Please contact the customer if the advanced user should be enabled, and ask the customer/admin what should be used as the advanced user password.
Also ask the customer/admin what should be used for other users. Then hand over the passwords to the customer/admin (see (→ New Password Handling in syngo / Page 94) and (→ Advanced User handling / Page 98))

3.1.17 MRAWP (MRC) Setup Auto Login Function



By default the autologin function is enabled. If this needs to be changed, see (→ New Password Handling in syngo / Page 94)

3.1.18 MRAWP (MRC) Performing COBR Image backup



Also, the user administrator does not exist in VE11D. The new administrative account is syadmin.



Please check the external medium (e.g. USB disk) first with the virus scanner of the service PC. Ensure that the external medium used is free of viruses.

ACTION	MRWP
■ Create the COBR image backup of the system accordingly - Manual Backup of Partition C to USB Disk or - Manual Backup of System Images to USB Disk Note: Due to performance and time issues, it is recommended to use a USB disk for backup in place of DVDs. The COBR image tool must be started in the following way due to Embedded Control used in VE11D: - Log in as user syadmin - Invoke a command shell from the Start menu (Start ⇒ All Programs⇒ Accessories ⇒ Command Prompt) - Type in the command: <code>C:\MedCom\bin\StartDiskImageGeneration.exe</code> ↵ Follow the instruction displayed in the cmd tool window.	☐

3.1.19 MRAWP (MRC) Enable the Auto Login Function



The following steps are required only if auto login is not enabled.

ACTION	MRAWP
■ Select Options --> Service --> Local Service .	☐
■ Enter the valid service key and click OK .	☐
■ Select Configuration in the service platform.	☐
■ Select local host Users	☐
■ Mark the checkbox " Enable Auto login " to enable the auto login function.	☐
■ Click Save .	☐
■ Click Home to leave the configuration.	☐
■ Close the local service UI.	☐

3.1.20 MRWP (MRSC) Startup



If you have a satellite console MRWP (MRSC) with shared database the frontend user (meduser) and administrative account (syadmin) and the passwords must be the same on both consoles MRAWP (MRC) and MRWP (MRSC).

Is this note the case the accounts frontend user (meduser) or/and administrative account (syadmin) at the MRAWP are locked after the reboot.

By by factory default the frontend user (meduser) and administrative account (syadmin) and the passwords are the same at both consoles.



If you want to change the frontend user name (meduser) or/and administrative account name (syadmin) or the passwords it is recommended to remove the network cable from the MRWP (MRSC). See also following SKB Entry (→ SKB0075323, <https://skb.healthcare.siemens.com/Pages/SKBAdvancedViewKBDocument.aspx?docid=SKB0075323>)



The Account lockout threshold defined in VE11C is the same as in VE11D. So the accounts for meduser and the administrative account (e.g. syadmin) are locked after using a wrong password four times. See also (→ New Password Handling in syngo / Page 94)



In some cases the following popups might be displayed. They can be ignored.

Fig. 60: Reboot popup


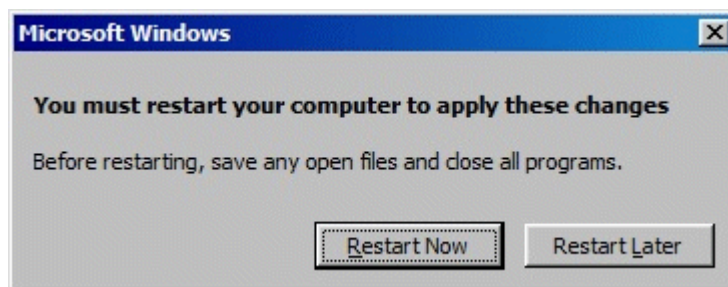


Fig. 61: Network location popup



ACTION	MRWP
▪ Switch on the MRWP.	☐
▪ If auto login is disabled, log in as meduser.	☐
▪ Wait until the main UI is displayed.	☐

3.1.21 MRWP (MRSC) Checking the Patient Browser Connection with Shared Database

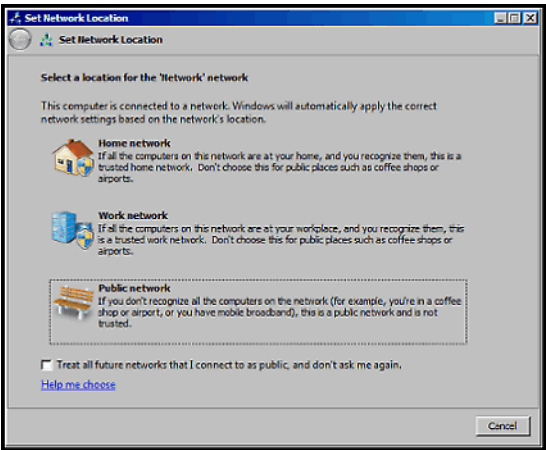

Ensure that the satcon with remote database configuration is not configured as a network node to the Main console. Otherwise problems will occur.

ACTION	MRWP
▪ Use the Patient Browser to check whether the MRAWP database can be reached.	☐

3.1.22 MRWP (MRSC) Network Printer Installation on the MRSC

ACTION	MRWP
▪ Install the network paper printer according to the software installation manual (this step is not part of the dynamic data) if previously installed on the MRAWP or MRWP.	☐
▪ If applicable, set the installed printer as the default printer. Select Start -> Settings -> Printers and Faxes Select the printer and click the right mouse button. In the menu displayed, select the "Default printer" entry.	☐

3.1.23 MRWP (MRSC) Configuration checks




With VE11D a new secure Service Key is introduced. More information can be found here (→ Service key handling with new Secure Service Key / Page 89)




The service key with copy and from the “*ServiceKey.dat” file. The file can be found in the following subdirectory C:\temp\license. Please delete this directory permanent after the update and store the service key file on a different location (e.g on the USB hard-disk which is used for the COBR image).



By default the advanced user is disabled. For the configuration check the advanced user must be enabled (see (→ MRWP (MRSC) Setup Advanced User Login Function / Page 80))

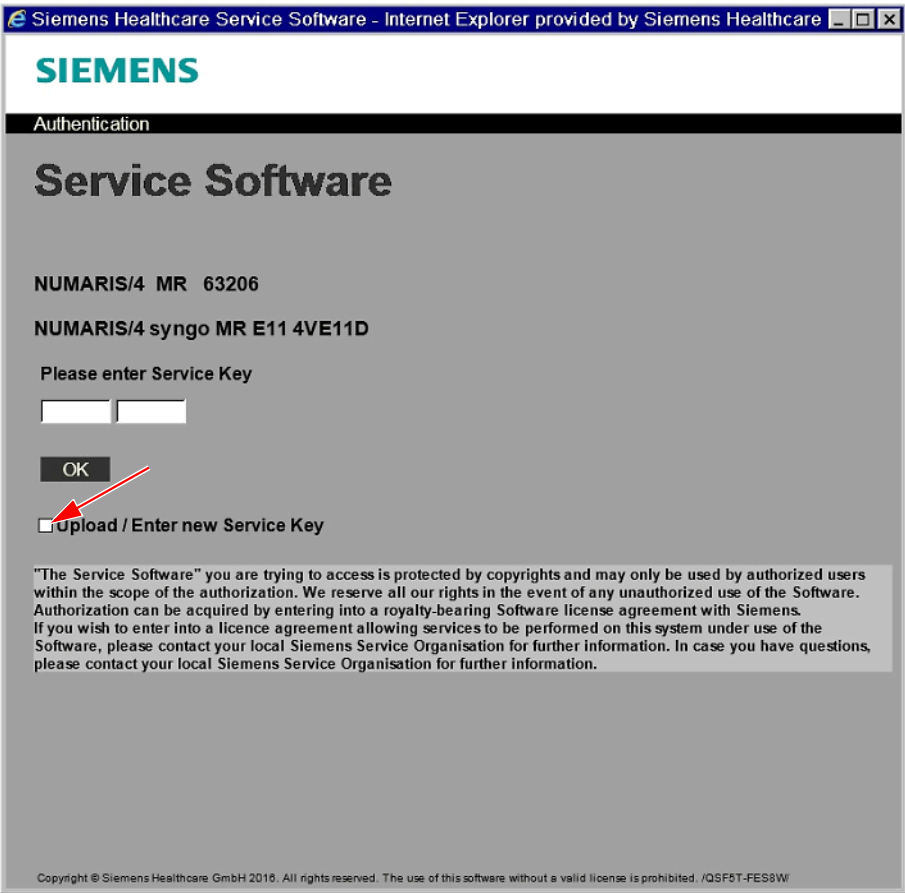ACTION	MRWP
Select Options ⇒ Service ⇒ Local Service.	☐
Select the checkbox “Upload/Enter new Service Key” to enable the upload <i>Fig. 62: Mark checkbox Upload</i> 	☐

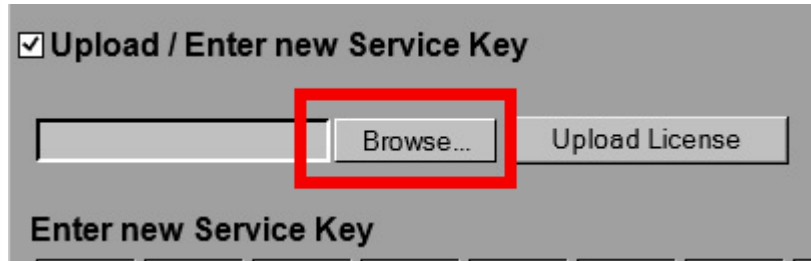ACTION	MRWP
Select the button Browse and navigate to the location of the "<code>*ServiceKey.dat</code>" file and select the file. <i>Fig. 63: Browse button</i>  The screenshot shows a dialog box titled 'Upload / Enter new Service Key'. It has a checked checkbox and a text input field. To the right of the input field are two buttons: 'Browse...' and 'Upload License'. The 'Browse...' button is highlighted with a red rectangular border.	☐
Select the button Upload License.	☐
Now the service key is uploaded and a login to the service UI is done.	☐
Select Configuration in the service platform.	☐


ACTION	MRWP
■ Select Monitor Type under Local Host and check the settings as follows: - Monitor Type: Color and Siemens (True Color) - Multi-Monitor settings Only "Start syngo in single monitor" is marked Fig. 64: Monitor settings SESO 	☐
■ Click Save if changes have been made and Finish	☐

ACTION

MRWP

- Under DICOM, select **General** and check if the checkbox Unencrypted is marked.

Fig. 65: Example from SIM Console

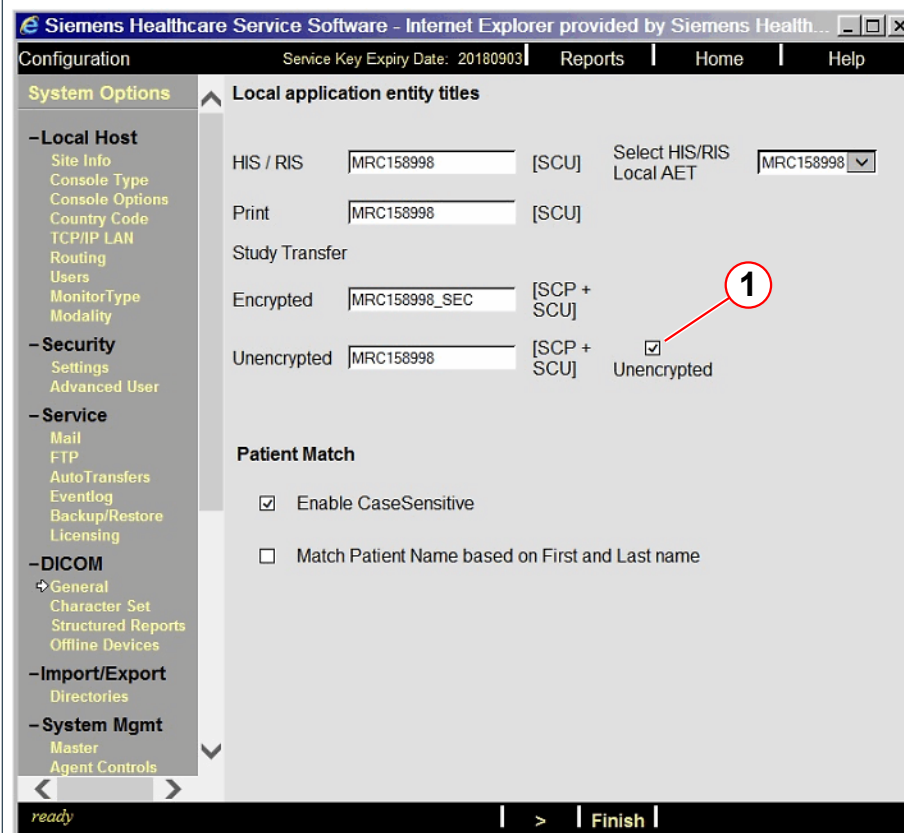

(1) Unencrypted checkbox must be marked

- Click **Save** if changes have been made and **Finish**

- Under DICOM, select **Network Nodes** and press under Host Properties >.

- On the next page, under **General Node Properties** in the drop-down menu **Logical Name**, select each available entry. For each entry check the **Connection Mode**. This must be set to **Unencrypted**.

Fig. 66: DICOM connection mode



Click **Save** if changes have been made for each entry.

All configured DICOM nodes from the former software version can be found under

C:\MedCom\MriSiteData\report.pdf or

C:\MedCom\MriSiteData\BaRe_HW_Info\MRC\report.html

If there are some missing, please carry out the configuration again.

- If all entries have been checked, click **Finish**.

ACTION	MRWP
Under DICOM, select Print Devices and press under Host Properties >.	☐
On the next page, under General Node Properties in the drop-down menu Logical Name, select each available entry. For each entry check the Connection Mode. This must be set to Unencrypted. <i>Fig. 67: DICOM connection mode</i>  Click Save if changes have been made for each entry.	☐
If all entries have been checked, click Finish.	☐
Click Home and confirm the requested restart of application with OK.	☐
Wait until the restart is performed and the syngo main UI is up and running.	☐



If a new Network Node, Print device or HIS/RIS Node is configured, be aware that the default setting for the Connection Mode is secure. Please change it to Unsecure.

3.1.24 MRWP (MRSC) Function Test of the DICOM Nodes

ACTION	MRWP
Send an image from the Patient Browser to a number of configured DICOM nodes. The images must be sent successfully	☐
Send a film sheet to the connected DICOM camera. The film sheet must be sent successfully to the DICOM camera.	☐



A DICOM connection can now be made secure (encrypted). This is a new configuration parameter. See (→ New feature for DICOM node configuration / Page 100)




The new DICOM Node Configuration Option "Supports Presentation State" must not be enabled for any DICOM storage node. The NUMARIS/4 VE11 does not currently support the "Presentation State" for created MR images.


3.1.25 MRWP (MRSC) System Management




Due to the CSSI process, system management cannot be activated at the factory. System management has to be activated at the customer site.

3.1.25.1 Installation and Activation Steps

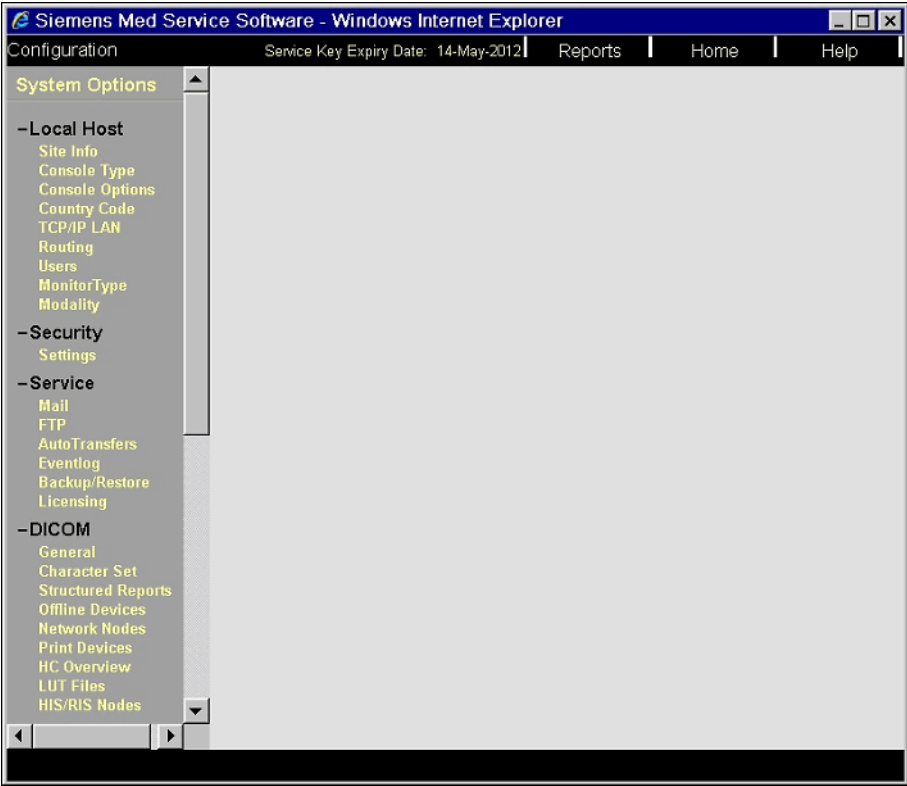ACTION	MRWP
Select Options --> Service -> Local Service.	☐
Enter the valid service key.	☐
In the local service UI, select Configuration	☐
Check that the settings for service region (Europe, USA, or Asia) under Local Host Site Info are correct. Correct the settings, if necessary, and save the changes. This entry assigns the IP addresses for the RDIAG server and the System Management server.	☐
Under the "Local Host" Routing, type in the IP address of the router or gateway used for SRS. In some cases, this entry is the default gateway that is used for the TCP/IP LAN configuration. See log file <code>C:\MedCom\MriSiteData\ExtBaRe\routing.txt</code> for the routing entry. For the installation of System Management, the correct entry under routing is essential. If this routing entry is not set (empty), you will get an error popup when entering the System Management configuration.	☐
Select System Options. Fig. 68: Configuration 	☐

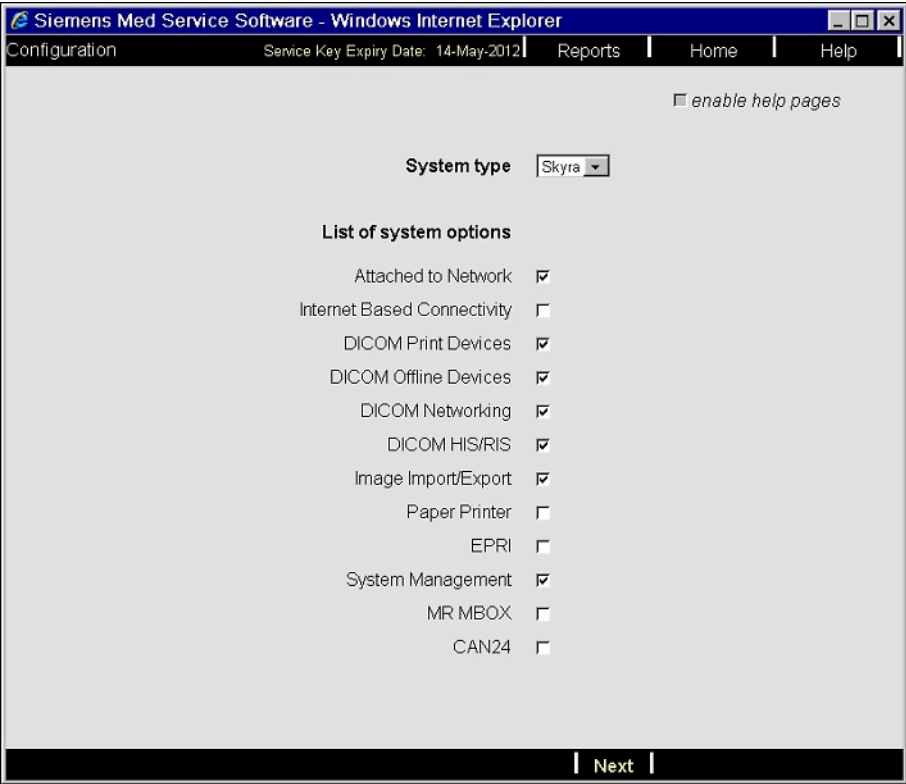ACTION	MRWP
■ Under "List of System Options", click the System Management checkbox <i>Fig. 69: System options</i> 	☐
■ Select Next.	☐

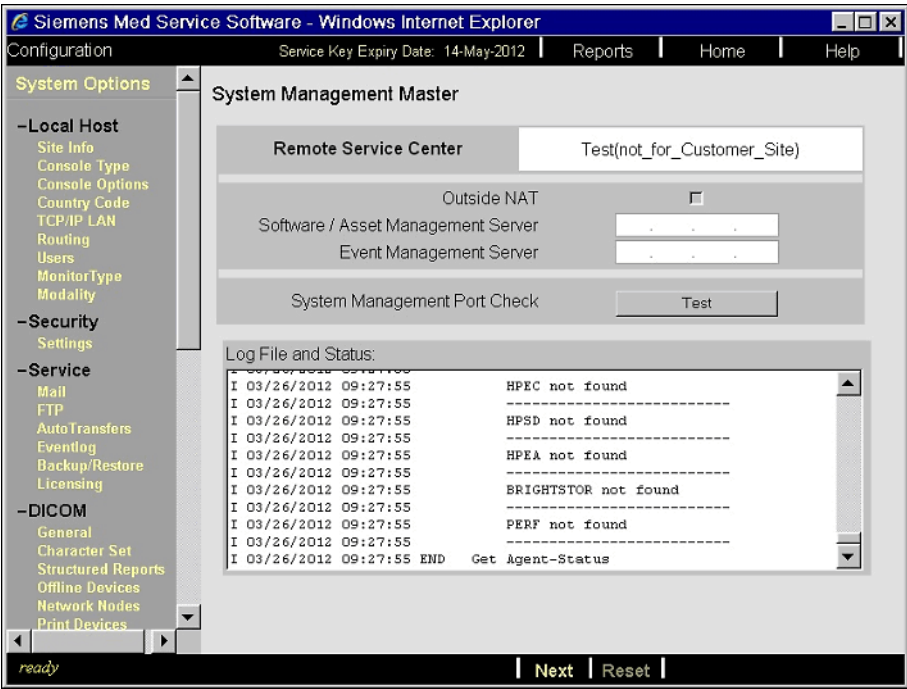ACTION	MRWP
■ On the System Management Master configuration page Fig. 70: System Management master page  See the report file C:\MedCom\MriSiteData\BaRe_HW_Info\MRC\report.html for the Outside NAT configuration.	☐


ACTION

MRWP

- Use the **Test** button to check whether the System Management server can be accessed.

☐

Fig. 71: System Management master can be reached

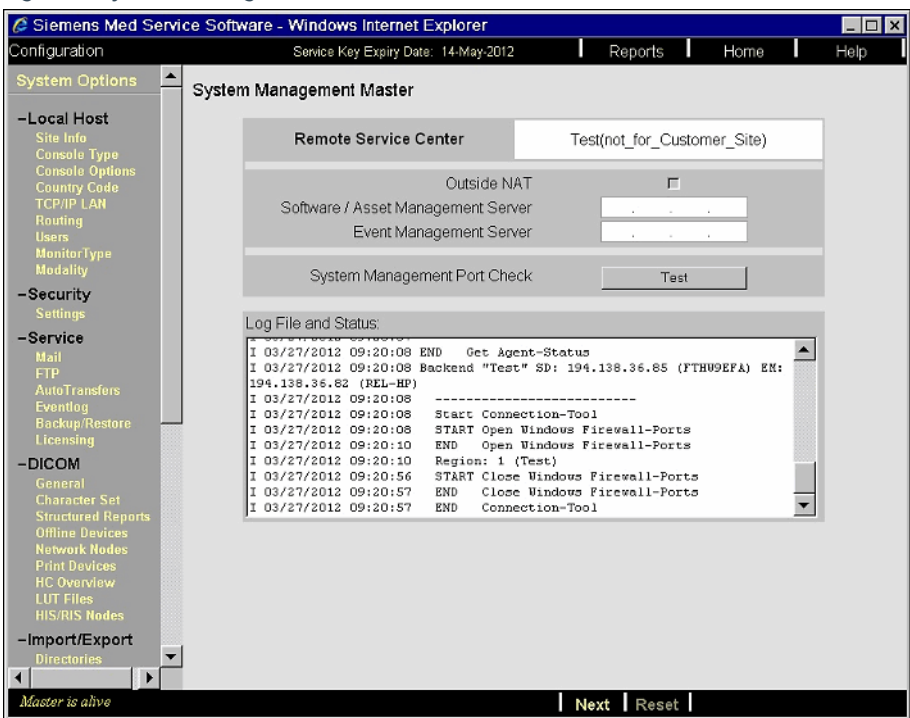

The log file for the test can be found under the directory:

C:\MedCom\service\mwtools\CCT\Logs

In the event of problems with the connection test, contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/ Technical Support Group.

- Select **Next**.

☐

ACTION

MRWP

- Under **System Mgmt. Agent Controls**, set the following configuration:
Mark the checkboxes "Go Active" and "Go Online".
(The other checkboxes have already been marked).

Fig. 72: System Management activation

The screenshot shows the 'System Management Functions' configuration window. It includes a table with columns 'Functions', 'Installed', and 'Status'. The 'Go Active' and 'Go Online' checkboxes are highlighted with a red circle and the number 1. Below the table, there is a 'Log File and Status' section showing a log of system management activities.

Functions	Installed	Status
SD- Asset- Management Agent	☒	not running
Event Management Agent	☒	not running
Go Active	☐	
Go Online	☐	

Log File and Status:

```
I 03/27/2012 08:07:32 -----  
I 03/27/2012 08:07:32 BRIGHTSTOR not found  
I 03/27/2012 08:07:32 -----  
I 03/27/2012 08:07:32 PERF not found  
I 03/27/2012 08:07:32 -----  
I 03/27/2012 08:07:33 END Get Agent-Status  
I 03/27/2012 08:07:33 Successful installation of System Management  
Base  
I 03/27/2012 08:07:33 -----  
I 03/27/2012 08:07:33 START Agent-Handling  
I 03/27/2012 08:07:33 END Agent-Handling
```

The System Management installation / registration / uninstallation may take several minutes. Please be patient while System Management is being installed / registered / uninstalled

(1) Mark both checkboxes

In the Log File and Status area, the current status of System Management activation/ installation is displayed.

- Click **Go**.

Note: The installation/registration will take approximately 10 to 15 minutes.

ACTION	MRWP
■ If activation is successful, the following is displayed: Fig. 73: System Management activation OK Select Finish to continue	☐
■ Click Finish and close the local service UI. In the event of problems (installation/activation fails a second time), contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/ Technical Support Group.	☐



For problems installing the Event Management Agent, see (→ System management installation workaround / Page 139)

3.1.26 MRWP (MRSC) Setup Auto Login Function



By default the autologin function is enabled. If this needs to be changed, see (→ New Password Handling in syngo / Page 94)

3.1.27 MRWP (MRSC) Setup Advanced User Login Function



By default the Advanced User Login function is disabled. If it needs to be changed, see (→ Advanced User handling / Page 98)



Please contact the customer if the advanced user should be enabled, and ask the customer/admin what should be used as the advanced user password.

Also ask the customer/admin what should be used for other users. Then hand over the passwords to the customer/admin (see (→ New Password Handling in syngo / Page 94) and (→ Advanced User handling / Page 98))

3.1.28 MRWP (MRSC) Perform a COBR Image Backup on the MRWP (MRSC)



Please check the external medium (e.g. USB disk) first with the virus scanner of the service PC. Ensure that the external medium used is free of viruses.

ACTION	MRWP
- Create the COBR image backup of the system accordingly - Manual Backup of System Images to USB Disk or - Manual Backup of System Images to DVD Note: Due to performance and time issues, it is recommended to use a USB disk for backup in place of DVDs. The COBR image tool must be started in the following way due to Embedded Control used in VE11D: - Log in as user syadmin - Invoke a command shell from the Start menu (Start ⇒ All Programs⇒ Accessories ⇒ Command Prompt) - Type in the command: <code>C:\MedCom\bin\StartDiskImageGeneration.exe</code> ↵ Follow the instruction displayed in the cmd tool window.	☐

3.1.29 MRWP (MRSC) Enable the Auto Login Function



The following steps are required only if auto login is not enabled.

ACTION	MRWP
Select Options --> Service --> Local Service.	☐

ACTION	MRWP
▪ Enter the valid service key and click OK .	☐
▪ Select Configuration in the service platform.	☐
▪ Select local host Users	☐
▪ Mark the checkbox " Enable Auto login " to enable the auto login function.	☐
▪ Click Save .	☐
▪ Click Home to leave the configuration.	☐
▪ Close the local service UI.	☐

3.1.30 Destroying Data on Old Host PCs



The data on old host PCs must be destroyed to make them unusable.
The old dynamic data sets (backup CDs) must also be destroyed.

3.1.30.1 BLANCCO, Erasure Software

Introduction

Warning

Before you start using the Blancco software, make sure that all files, folders, software applications or any other information that you want to save for later use are backed up on an appropriate medium other than the original hard drive. Data that has been erased from the hard disk by this program cannot be recovered by any existing method.



WARNING

Selecting the wrong hard disk can result in
loss of data with the use of blancco.

» Ensure that the hard disk selected for deletion is the correct one.

Blancco Information

Blancco Server Edition can erase up to 16 hard disks simultaneously

Creating the blancco CD

An ISO image of the **blancco** CD is available for download from: Product Information ⇒ Common Services ⇒ Service Tools

Userdocumentaion (→ Blancco Manual, <https://productinformation.healthcare.siemens.com/FolderFile.ashx?filePath=XFNlcnZpY2UgVG9vbHNcQmxhbmNjb19Ecml2ZV9FcmFzZXJfTWfudWFsX2VuLVVTLnBkZg%3d%3d&uuid=eda948cf-7154-428c-a3ca-ef569f6a0ec7>)

Iso file (→ Blancco Drive Eraser, <https://productinformation.healthcare.siemens.com/FolderFile.ashx?filePath=XFNlcnZpY2UgVG9vbHNcRHJpdmVFcmFzZXJfc2VydmVyXzYxMF9SU192MzEwMzE4Lmlzbw%3d%3d&uuid=eda948cf-7154-428c-a3ca-ef569f6a0ec7>)

1. Download the ISO image.
2. Label the blank CD as "**blancco** Erasure Software"
3. Insert the CD into the drive of your service computer.
4. Use a standard burn program (such as Nero) to open the file and burn the CD image to the CD.
5. Remove the CD from the drive.

Using the blancco CD


Bootting with blancco

1. Connect monitor, mouse, keyboard, and power cord.
2. Switch on the PC.
3. Insert the "**blancco** Erasure Software" CD into the PC with the drive(s) to be deleted.
4. Switch off the PC.
5. Switch on the computer power to boot the system from the CD (if necessary, change the BIOS boot priority - to do this you have to access the BIOS during the boot process by pressing the **F2** or **F8** key).
 - » The **blancco** home screen is visible with all identified hard disks of the PC.



If necessary, press **F2** to change the language and the keyboard layout.

Fig. 74: Blanco language / keyboard layout

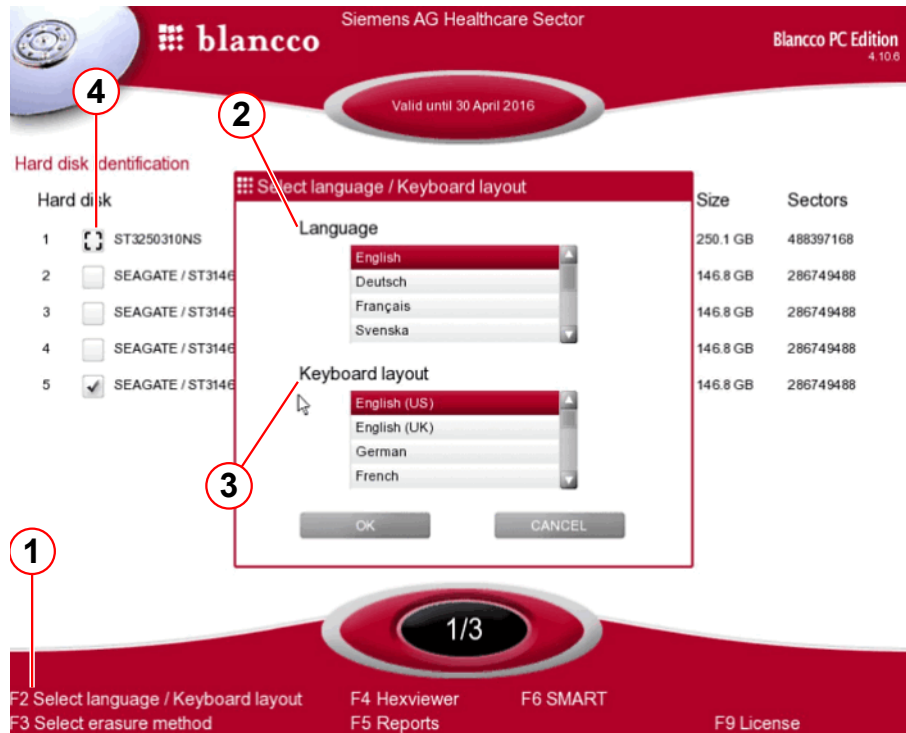

- (1) F2 to open the menu
- (2) Language
- (3) Keyboard layout
- (4) Selection of hard disks

Deletion with blanco


1. Select the hard disks to be deleted (→ 4/ Fig. 74 Page 83).



blanco can erase up to 16 hard disks simultaneously.

2. Press **Erase**.
 - » A warning message appears.

Fig. 75: Erase disk(s) warning

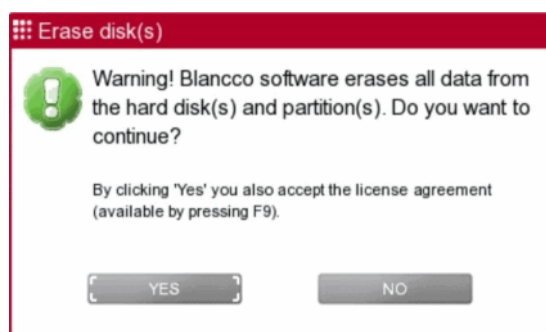

3. Press **Yes** to start the erasure process.
 - » The erasure process starts and the progress and duration are visible.
4. Click on **OK** when the process is finished and the **Erase successful** message appears.
5. Insert a memory stick to save the report.
6. Press **Save** or **F10**.
 - » The report is saved to the memory stick.
7. Press **Ok**.
 - » The process is finished and the PC can be turned off.

3.1.31 Final Check - to be performed for MRAWP (MRC) and MRWP (MRSC) (if applicable)



If the customer uses a Network Time Protocol Server for synchronizing the time on the MR system, the configuration has to be done again. This configuration cannot be done in the factory due to the check to determine if the NTP server is reachable (local service UI, Site Info, step 2)

- Verify the correct installation as follows:
In the main syngo MR UI, select
Help -> Info
Check whether the following entry is visible under the Siemens software packages
- N4_VE11D_LATEST_20180608
The label
- N4_VE11D_LATEST_20180608
indicates a correct VE11C base installation.
- Check the RDIAG configurations on the **MRAWP (MRC)** and, if applicable, on the **MRWP (MRSC)** as well (the MRWP (MRSC) must also be connected to RDIAG) as follows:
 - Check the host name and IP address of the configured RDIAG server (Local Service -> Configuration -> Service Mail)

	Host name	IP address
Europe, Africa	LUX09505	194.138.39.18

	Host name	IP address
North/South America	SRSACC1	129.73.116.92
Asia	sgpt806x	194.138.243.178

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- In the second step of the mail configuration, check the defined mail target (for the STE, the e-mail address must be changed according to the following table)

	Logical name for FTP target	Select the mail server host	E-mail address
Europe, Africa	SRS_mail	LUX09505	mrrsmc.team@siemens-healthineers.com
North/South America	SRS_mail	SRSACC1	mrrsmc.team@siemens-healthineers.com
Asia	SRS_mail	sgpt806x	mrrsmc.team@siemens-healthineers.com

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- In the second step of the FTP configuration (Local Service -> Configuration -> Service FTP), check the defined FTP target

	Logical name for FTP target	FTP host selection
Europe, Africa	SRS_ftp	LUX09505
North/South America	SRS_ftp	SRSACC1
Asia	SRS_ftp	sgpt806x

FTP user name: `rdiaftp`

Password: `siemens`

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Test the FTP connection to the RDIAG server as follows:

Local Service -> File&Image Tools

Source: `File`

Action: `Transfer to System`

File Type: `Any file`

From source: `SRS_ftp`

If a list from the directory `/outgoing/mr` (directory on the RDIAG server) is displayed, the connection and FTP access are OK

If the test is not OK, check the SRS configuration; refer to

- Click Home

- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **EL_Sample**. The following settings must be applied:
Active checkbox: must be marked
Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).
- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **EL_Transfer**. The following settings must be applied:
Active checkbox: must be marked
Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).
- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **IB_Sample**. The following settings must be applied:
The "Activated" checkbox must be marked
Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).
- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **IB_Transfer**. The following settings must be applied:
The "Activated" checkbox must be marked
Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).
- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **edb**. The following settings must be applied:
Active checkbox: must be marked
Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).
- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **qdb**. The following settings must be applied:
Active checkbox: must be marked
Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).
- The **UTILDATA** auto report is used for transferring the utilization data of the system. Utilization Management is optional; for legal reasons the customer has to sign a contract in advance. Please consult your service organization before you activate System Utilization.



Prevent failure messages where the regional language setting did not match, e.g. UI language

- guidance text used within the workflow (Siemens pre-defined texts)
- language text used in error messages
- table position of the protocol after opening in a completely different range, and not useful for planning. (Due to mismatch between "." and ";")

See SKB0022441 entry and (→ Failure messages where the regional language setting did not match, e.g. UI language / Page 105).

3.1.32 Customer Information

The mandatory (5 days on site) Appl. Training is being provided via MR00x/16/E.
Any additional training required has to be determined individually.

3.1.33 Removed Windows Feature in VE11D

Due to OS hardening, the following Windows features are removed from the VE11D.

- Windows Gadget platform
- Windows search
 - Search functionality in task menu
 - Search functionality in Windows Explorer
- XPS Viewer
- Media features
 - Windows DVD maker
 - Windows Media Center
 - Windows Media Player
- Remote Differential Compression
- Tablet PC Components
- Games
- Simple TCP Services
- TFTP Client

3.1.34 Final Work Steps

1. **For systems equipped with the Expert-i license:**
Follow the instructions in document MR00-000.814.28 to update existing client PCs using the Expert-i client software delivered with the CSSI upgrade.
 Otherwise scrap the Expert-i client software.
2. Update the system documentation. Fill out, and if necessary, make a copy of the attached "Completion Protocol" and file it in the corresponding folder.

3. The update has to be reported by authorized personnel, either via an intranet application or by automated data transfer (only for countries connected via the update handling interface).
4. **Return of parts:** Please return the hardware that was replaced by the HW/SW Upgrade if the customer agrees that Siemens should take care of the disposal of the old HW.
5. **IVK Update:** New IVKs installed as a result of an update are to be assigned to the system-identifying IVK.

As required, replaced IVKs are to be assigned to user status SCRP (scrapped) or RETG (returned goods), depending on the agreement with the customer. Please clarify with your local country data management if necessary.

3.1.35 Changes to Previous Version

n.a.

4.1 Appendix

4.1.1 Service key handling with new Secure Service Key



The new Service Key is now 280 characters long (56 x 5 characters without dash)
The service software login page offers by default the input field for the last two character groups (last 10 characters). The rest of the key is stored on the system during the initial key upload. Every new key must be uploaded first. The system evaluates the service level automatically from the entered last 10 characters.

4.1.1.1 syngo Service Software Login with service key import

Service Key Upload from file



The service key upload from the “*ServiceKey.dat” file is only possible, if not more than one service key level is in the “*ServiceKey.dat” file. Are more service key levels in the file the upload may fail.

- Select the checkbox **"Upload/Enter new Service Key"** to enable the upload

Fig. 76: Mark checkbox Upload

Siemens Med Service Software - Internet Explorer

SIEMENS

Authentication

Service Software

PN SW 0

syngo VH22C SL24P29 VE10A

Please enter Service Key

Enter Serial No

OK

☐ Upload / Enter new Service Key

"The Service Software" you are trying to access is protected by copyrights and may only be used by authorized users within the scope of the authorization. We reserve all our rights in the event of any unauthorized use of the Software. Authorization can be acquired by entering into a royalty-bearing Software license agreement with Siemens.
If you wish to enter into a licence agreement allowing services to be performed on this system under use of the Software, please contact your local Siemens Service Organisation for further information. In case you have questions, please contact your local Siemens Service Organisation for further information.

- If no configuration was done before, the the system serial number must be typed in the field **"Enter Serial No."**

Fig. 77: Enter Ser. number

Please enter Service Key

Enter Serial No

OK

☒ Upload / Enter new Service Key

- Select the button **Browse** and navigate to the location of the "`*ServiceKey.dat`" file and select the file.

Fig. 78: Browse service key location



- Select the button **Upload License**
- Now the service key is uploaded and a login to the service UI is done.

Service Key Upload with copy and paste



The service key with copy and from the "`*ServiceKey.dat`" file in every case possible.

- Select the checkbox **"Upload/Enter new Service Key"** to enable the upload

Fig. 79: Mark checkbox Upload

Siemens Med Service Software - Internet Explorer

SIEMENS

Authentication

Service Software

PN SW 0

syngo VH22C SL24P29 VE10A

Please enter Service Key

Enter Serial No

0

OK

☐ Upload / Enter new Service Key

"The Service Software" you are trying to access is protected by copyrights and may only be used by authorized users within the scope of the authorization. We reserve all our rights in the event of any unauthorized use of the Software. Authorization can be acquired by entering into a royalty-bearing Software license agreement with Siemens.
If you wish to enter into a licence agreement allowing services to be performed on this system under use of the Software, please contact your local Siemens Service Organisation for further information. In case you have questions, please contact your local Siemens Service Organisation for further information.

- If no configuration was done before, the the system serial number must be typed in the field **"Enter Serial No."**

Fig. 80: Enter Ser. number

Please enter Service Key

Enter Serial No

92116


OK

☒ Upload / Enter new Service Key

- Open the **"*ServiceKey.dat"** file with wordpad and go the the section with your service leve e.g. Servicelevel 7.

- Mark only the characters, which are marked in the figure below and press the keys Ctrl + C in order to copy the selection.

Fig. 81: Copy Paste Service Key

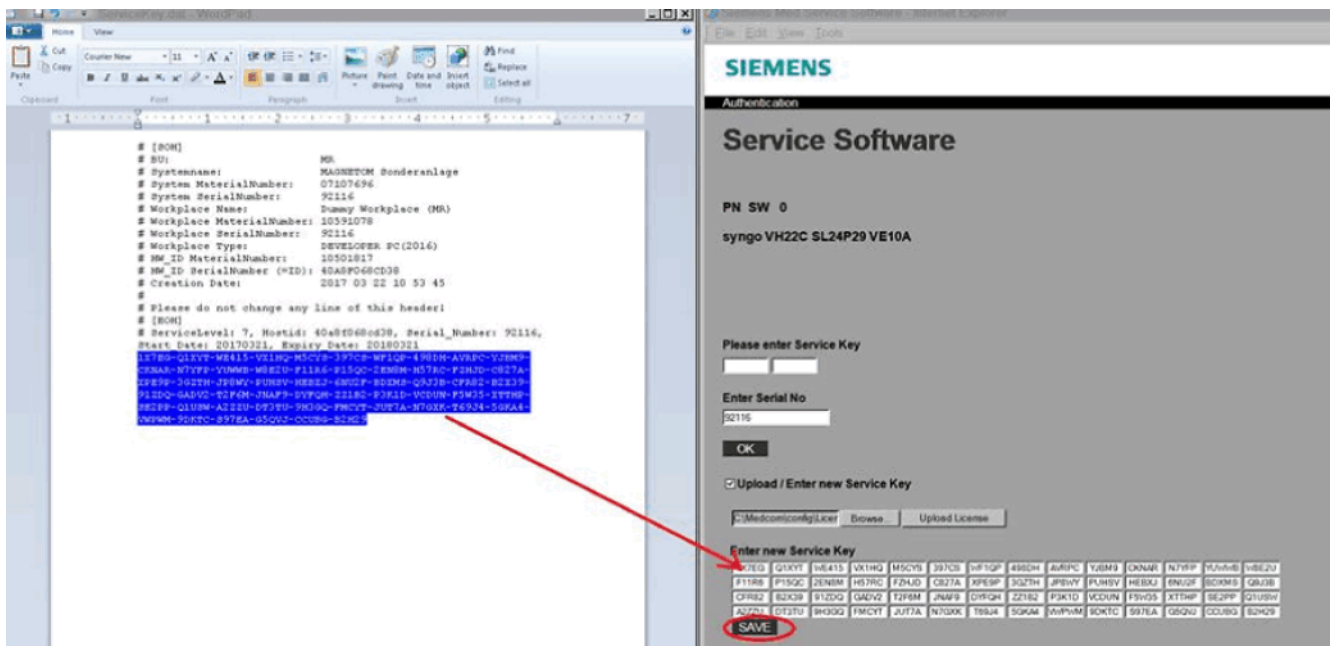

- Select the first field under “Enter new Service Key” and press the keys CTRL + V in order to paste the whole service key.
- Select **Save**.
- Now the service key is uploaded and a login to the service UI is done.


4.1.1.2 syngo Service Software Login



The service software login page offers by default the input field for the last two character groups (last 10 characters). The rest of the key is stored on the system during the initial key upload. Every new key must be uploaded first. The system evaluates the service level automatically from the entered last 10 characters.

- Select in the syngo main UI: **Options** ⇒ **Local Service**
- Enter the last 10 characters from the service key into two shown groups.

Fig. 82: Service Key groups

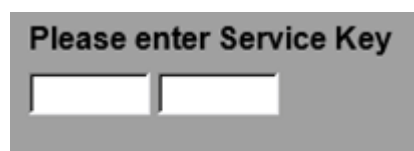

This can be done either by typing in each character or via copy/paste from the service key file.

- Press **OK** to login.

4.1.2 New Password Handling in syngo



There are no longer fixed or default passwords within syngo. Only random or individually defined passwords are valid. OS Hardening is enabled with Numaris VE11C. With OS Hardening enabled, a specific Administrator account has to be defined. The default "SYADMIN" is suggested. This account is created during the installation in the factory. Autologin for the meduser is also enabled during the factory installation. The default passwords are the same as in VE11C for the frontend user (meduser) and administrative account (syadmin).



The Account lockout threshold defined like in VE11C is the same in VE11D. So the accounts for meduser and the administrative account (e.g. syadmin) are locked after using a wrong password four times.

4.1.2.1 Disable Autologin

ACTION	MRAWP/MRWP
■ Select Options ⇒ Service ⇒ Local Service .	☐
■ Enter the valid service key (service level 4 or 7) and click OK .	☐
■ Select Configuration in the service platform.	☐
■ Select Local Host Users	☐
■ Remove the Enable Autologin checkbox mark to disable the Autologin function.	☐
Fig. 83: SESO users 	

ACTION	MRAWP/MRWP
- The password for both users (frontend user (e.g. meduser) and administrative account (e.g. syadmin)) must now be set. The password must meet the following conditions: - Minimum length of 14 characters - It should have 2 uppercase letters and 2 lowercase letters - It should have 2 numbers and 2 special characters.	☐
Click Save.	☐
Click Finish.	☐
Click Home to leave the configuration.	☐
Close the local service UI	☐




See also (→ Microsoft info, [https://technet.microsoft.com/en-us/library/cc786468\(v=ws.10\).aspx](https://technet.microsoft.com/en-us/library/cc786468(v=ws.10).aspx)) for password requirements and allowed characters.

4.1.2.2 Enable Autologin

- If on site with local access to the service software, and if password for the meduser account is known

ACTION	MRAWP/MRWP
Select Options ⇒ Service ⇒ Local Service.	☐
Enter the valid service key (service level 4 or 7) and click OK.	☐
Select Configuration in the service platform.	☐
Select Local Host Users	☐

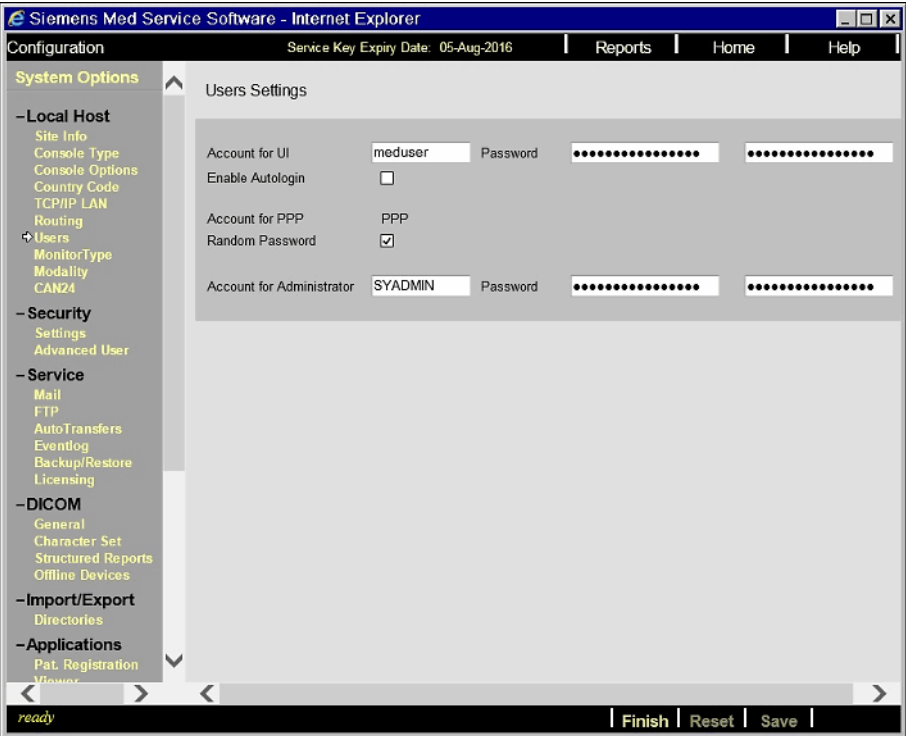ACTION	MRAWP/MRWP
Mark the Enable Autologin checkbox to enable the Autologin function. Fig. 84: SESO users 	☐
Click Save.	☐
Click Finish.	☐
Click Home to leave the configuration.	☐
Close the local service UI	☐

- If AutoLogin is disabled and meduser password is not known, use Service Software Utilities from Remote (Limited Permanent) to enable AutoLogin.
To do this, LAN or SRS connection must be possible (IP or host name is known, Service Key)

ACTION	MRAWP/MRWP
Remote login via SRS or LAN.	☐
Enter the valid service key (service level 3 or 7) and click OK.	☐
Select Utilities in the service platform.	☐

ACTION	MRAWP/MRWP
On the Utilities page, under Source, select the entry Enable Autologin. <i>Fig. 85: Source entries</i> 	☐
Click Go. <i>Fig. 86: Utilities Autologin success</i> 	☐
Click Home.	☐
Close the service UI	☐


Now autologin is enabled and the password can be set again.

4.1.3 Advanced User handling

4.1.3.1 Enable Advanced User



The login for the advanced user can now be enabled or disabled. Also a password must be set if the advanced user is enabled in the Service Software UI.

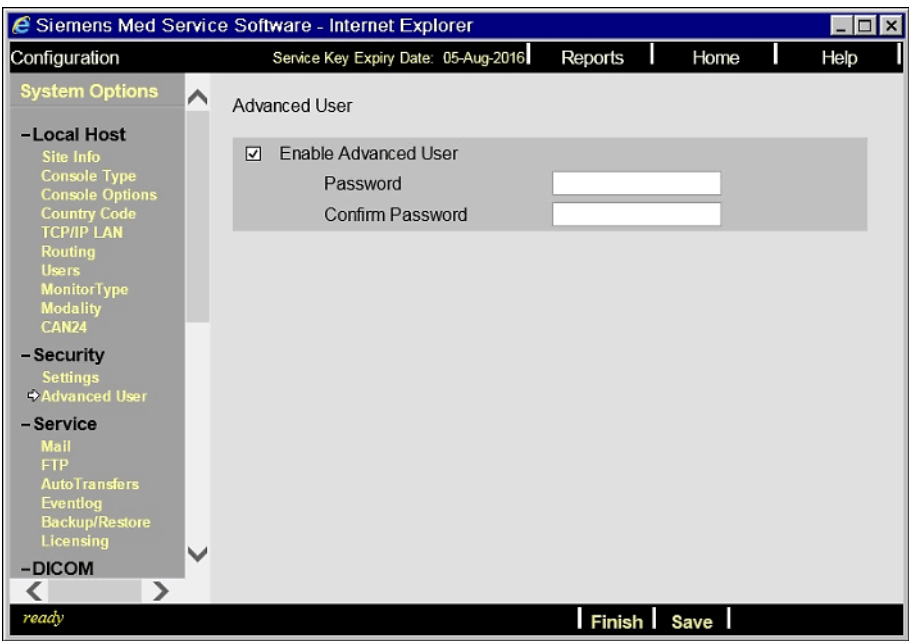ACTION	MRAWP/MRWP
■ Select Options ⇒ Service ⇒ Local Service .	☐
■ Enter the valid service key (service level 7) and click OK .	☐
■ Select Configuration in the service platform.	☐
■ Select Security Advanced User	☐
■ Mark the Enable Advanced User checkbox to enable the Advanced User login in the Windows Start menu. <i>Fig. 87: SESO Advanced User</i> 	☐
■ The password for the advanced user must now be set. The password must meet the following conditions: - Minimum length of 14 characters - It should have 2 uppercase letters and 2 lowercase letters - It should have 2 numbers and 2 special characters. The following 4 special characters are not allowed for the advanced user password: \ * ? ^	☐
■ Click Save .	☐
■ Click Finish .	☐

ACTION	MRAWP/MRWP
▪ Click Home to leave the configuration.	☐
▪ Close the local service UI	☐

4.1.3.2 Advanced User login

Procedure to display the complete Start menu:

- Run the Advanced User program from the Start menu (<CTRL> + <ESC>):
Start -> Programs -> Advanced User.
- Enter the user password:

You now have access to the complete Start menu.


4.1.3.3 Disable Advanced User



The login for the advanced user can now be enabled or disabled. No password must be set if the advanced user is disabled in the Service Software UI.

ACTION	MRAWP/MRWP
▪ Select Options ⇒ Service ⇒ Local Service .	☐
▪ Enter the valid service key (service level 7) and click OK .	☐
▪ Select Configuration in the service platform.	☐
▪ Select Security Advanced User	☐
▪ Remove the mark from the checkbox Enable Advanced User to disable the Advanced User login in the Windows Start menu.	☐

Fig. 88: SESO Advanced User

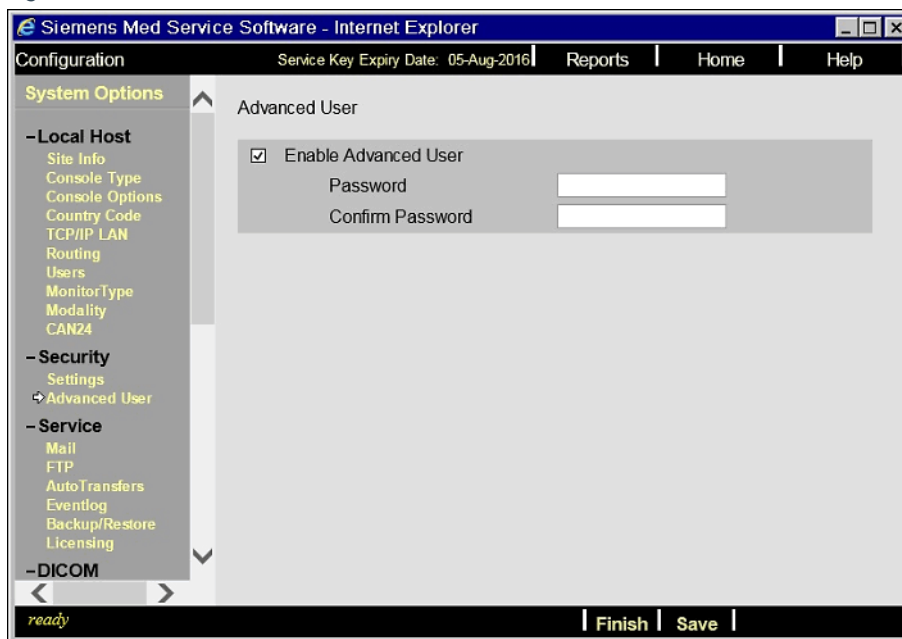


ACTION	MRAWP/MRWP
■ Click Save .	☐
■ Click Finish .	☐
■ Click Home to leave the configuration.	☐
■ Close the local service UI	☐

4.1.4 New feature for DICOM node configuration

With the new syngo version VH22B it is possible to configure the connection mode for each DICOM node (Network Nodes, Print Devices, HIS/RIS Nodes). If under connection mode the radio button secure (default) is selected, the DICOM communication is using the SSL/TSL protocol for the connection.

Fig. 89: DICOM Nodes configuration

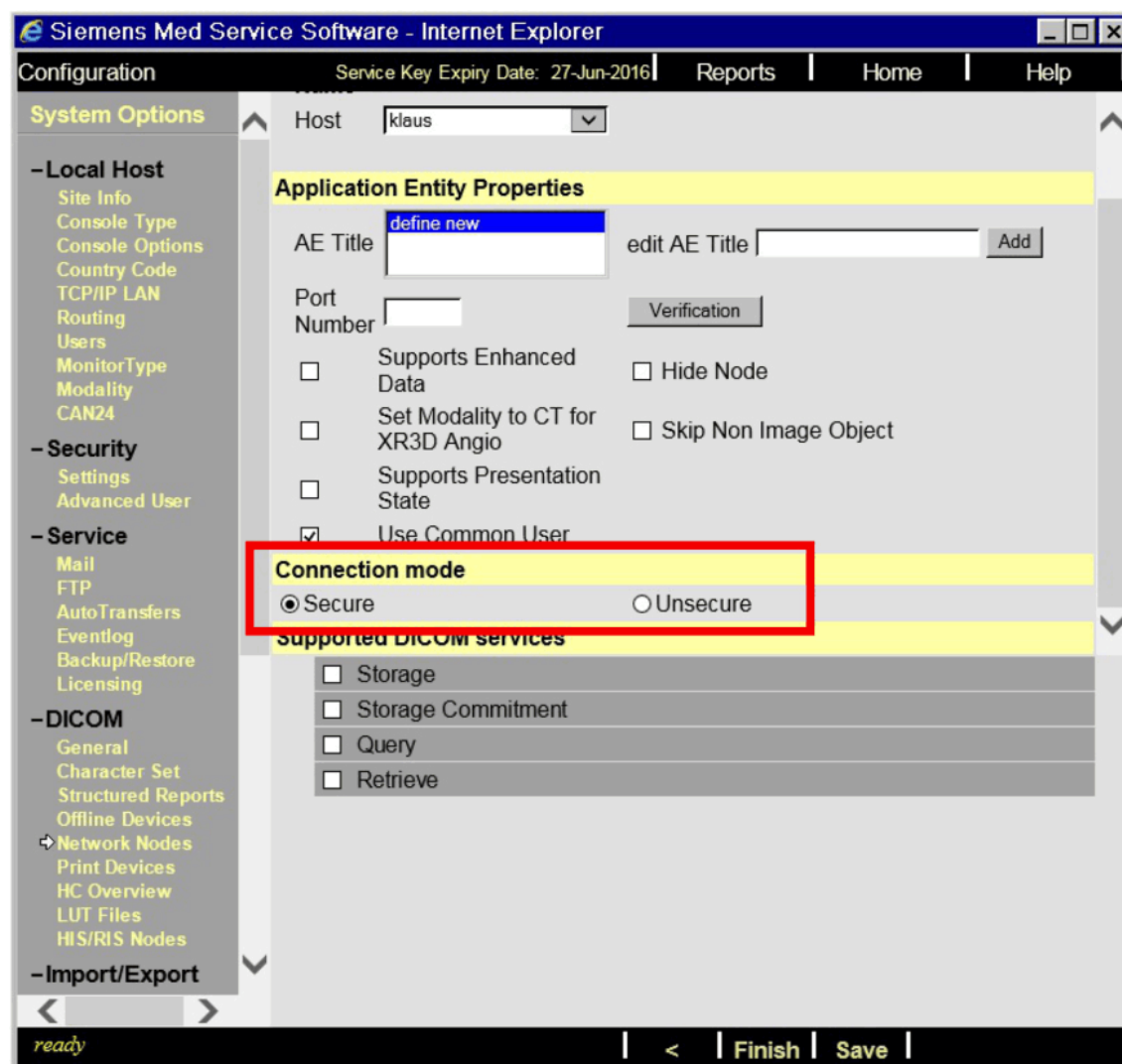

The connection fails if the DICOM node is not supporting the secure DICOM connection. In the event log the following message can be found (example).

151030|17:43:45|CSA_DSA|16|TLS Handshake started with CONNECTOR. Proposing Cipher suite AES128-SHA:DES-CBC3-SHA

151030|17:43:45|CSA_DSA|E17|TLS Handshake rejected by CONNECTOR.

151030|17:43:45|CSA_DAL|E1|Error in Dicom Association Layer: An error occurred in Application 1001 within Association 0

In this case modify the connection mode from secure to unsecure for this DICOM node.



Please be aware that the default for the connection mode is secure if a new DICOM node is configured in the Service UI under configuration.

Fig. 90: Default

Connection mode

☒ Secure

☐ Unsecure

The connection mode configuration must be changed to unsecure. This is currently the general use case for the DICOM node configuration.

Use secure mode only if the DICOM node supports the secure DICOM communication with the SSL/TSL protocol.

4.1.5 Drive Letter Assignment for the CD/DVD-ROM Drive and/or for the CD/DVD-R Drive

4.1.5.1 Description:

The drive letter assignment for the CD/DVD-ROM drive and/or for the CD/DVD-R drive is not correct. The following is displayed in the local service configuration

Fig. 91: Error message



Fig. 92: Error message


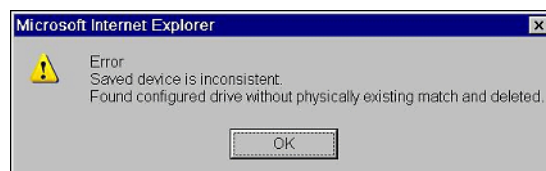


Fig. 93: Local service UI

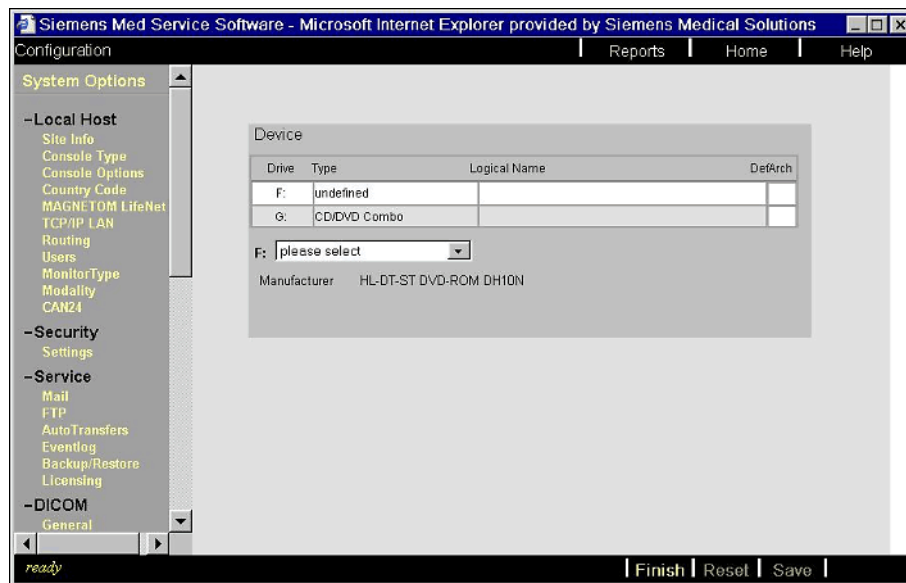


4.1.5.2 Solution:

The following procedure describes the correction for the drive letter assignment.

- Log in as the "syadmin".
- Start the disk manager by selecting **START** and entering `diskmgmt.msc` ↵


Fig. 94: Start Menu

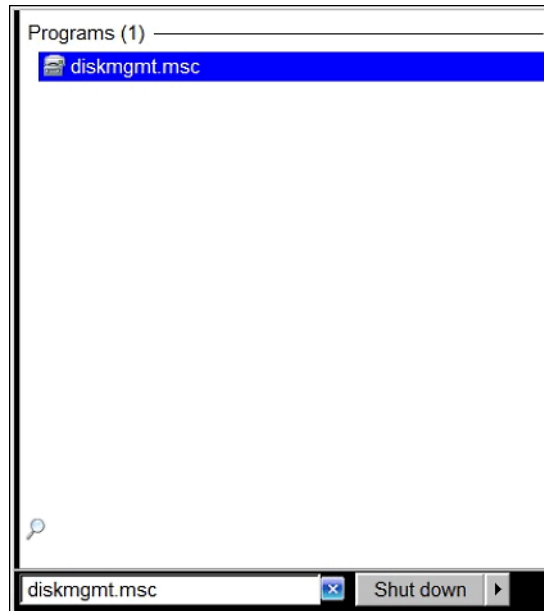

- The disk manager pops up.
- Select the drive to be (re-) assigned with a right mouse click, select **'Change drive letter and paths...'** from the popup menu, select **'Change'** in the next popup and choose the proper drive letter from the drop-down list of the next one. Click **'Yes'** in the 'Con-

firm' popup. (Repeat this sequence for every drive for which you have to change the drive letter).

Fig. 95: Disk management

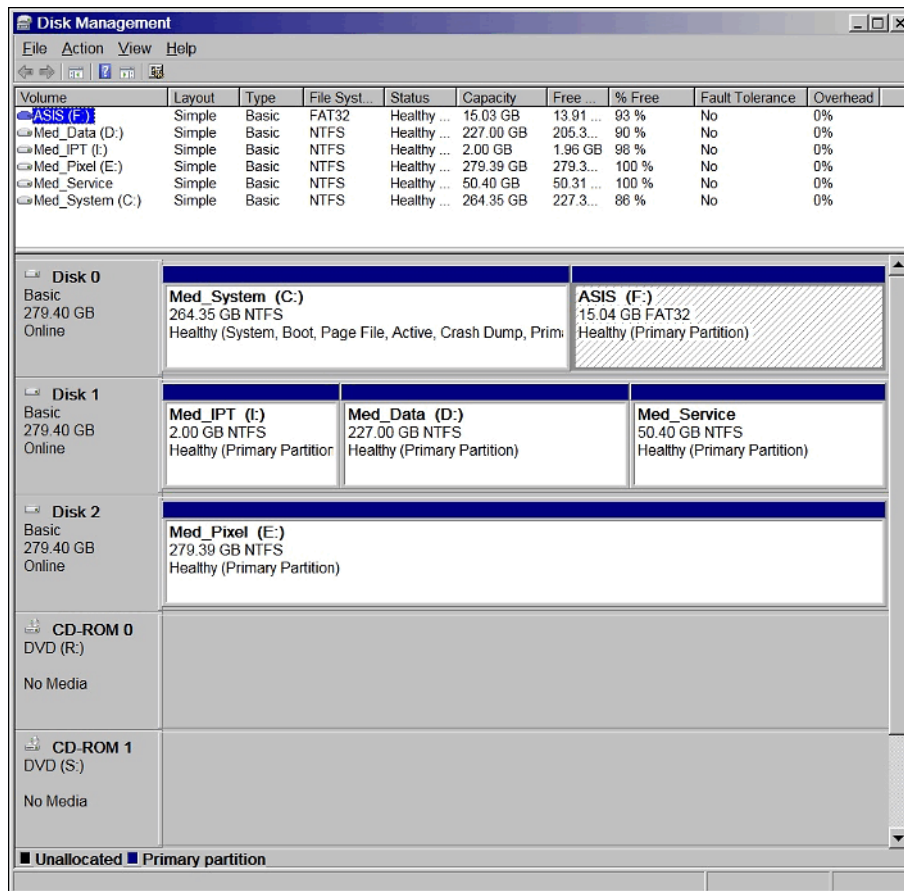

- CD/DVD-ROM drive -> R:
- CD/DVD-R (writer) drive-> S:
- Perform a reboot.
- Reconfigure Offline devices

4.1.6 Correct syngoUI Language Settings

Problem: No message displayed in the MR Scanner & Acquisition popup.


The reason is that the syngo main UI language is not correct set, wrong format is used. Please correct it see below.

Language change:

- Select in the syngo main UI: **Options** ⇒ **Configuration** ⇒ **Regional and Language Options**

Select the correct language in the Formats tab under the drop down box Format.

Fig. 96: Regional and Language Options

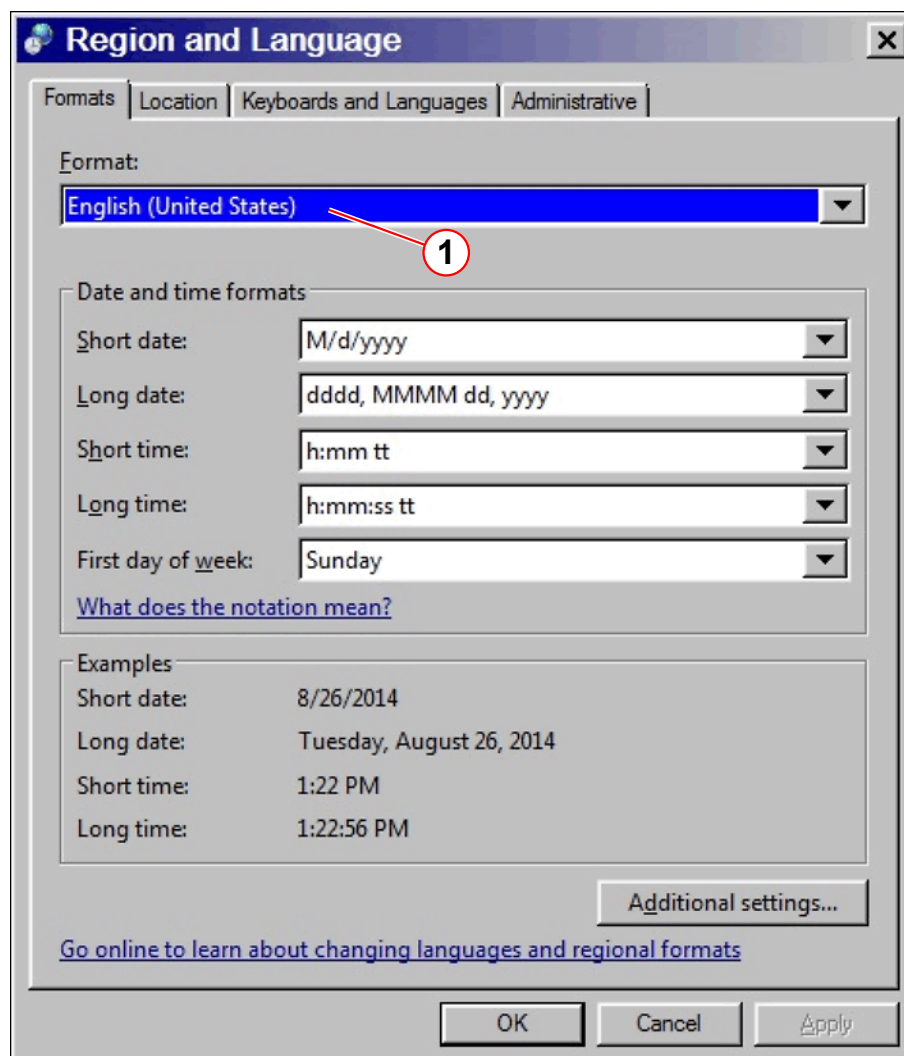

Only the following entries can be used for 1 in the graphic above:

Fig. 97:

main UI language	entry in the pull down menu
german	German (Germany)
english	English (US)
english (metric)	English (UK)
french	French (France)
spanish	Spanish (Spain)
chinese	Chinese (PRC)
japanese	Japanese

- Press OK for restart the system
- Check the scanner symbol. If the MARS isn't connected, the scanner symbol has to be crossed out. Select this symbol, and a new window with messages should appear

4.1.7 Failure messages where the regional language setting did not match, e.g. UI language

The following failure messages can occur:

- guidance text used within the workflow (Siemens pre-defined texts)
- language text used in error messages
- table position of the protocol after opening in a completely different range, and not useful for planning. (Due to mismatch between "." and ",")

The following has to be checked in order to prevent these mentioned issues:

- Before doing this: Close the patient!
- Deactivate Advanced User! (Start > All Programs > Startup > MR_User)

- Select from the syngo main UI Options > Configuration > Regional Settings > use the first **format** entry

Fig. 98: Example for English


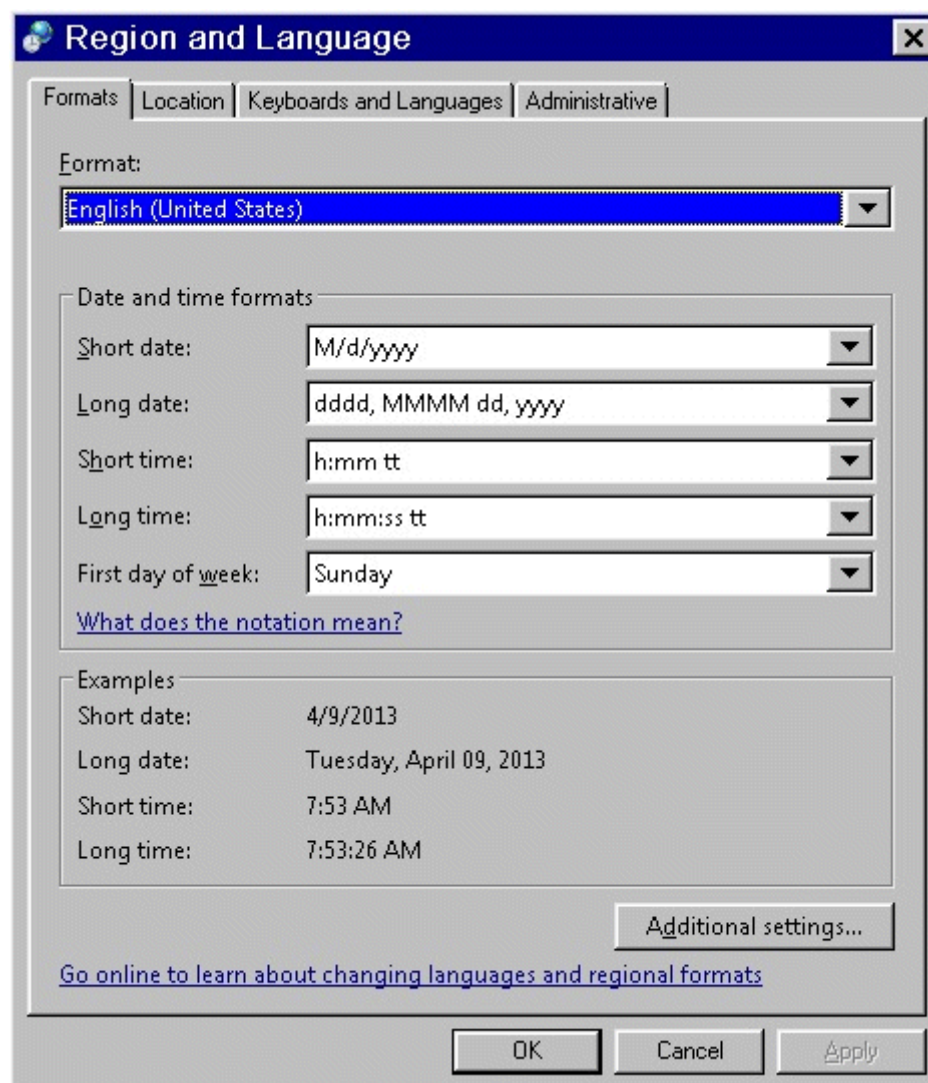


Fig. 99: Example for German



Press **Apply** and **OK**.

- After switching the language a reboot is pending. This reboot is mandatory!

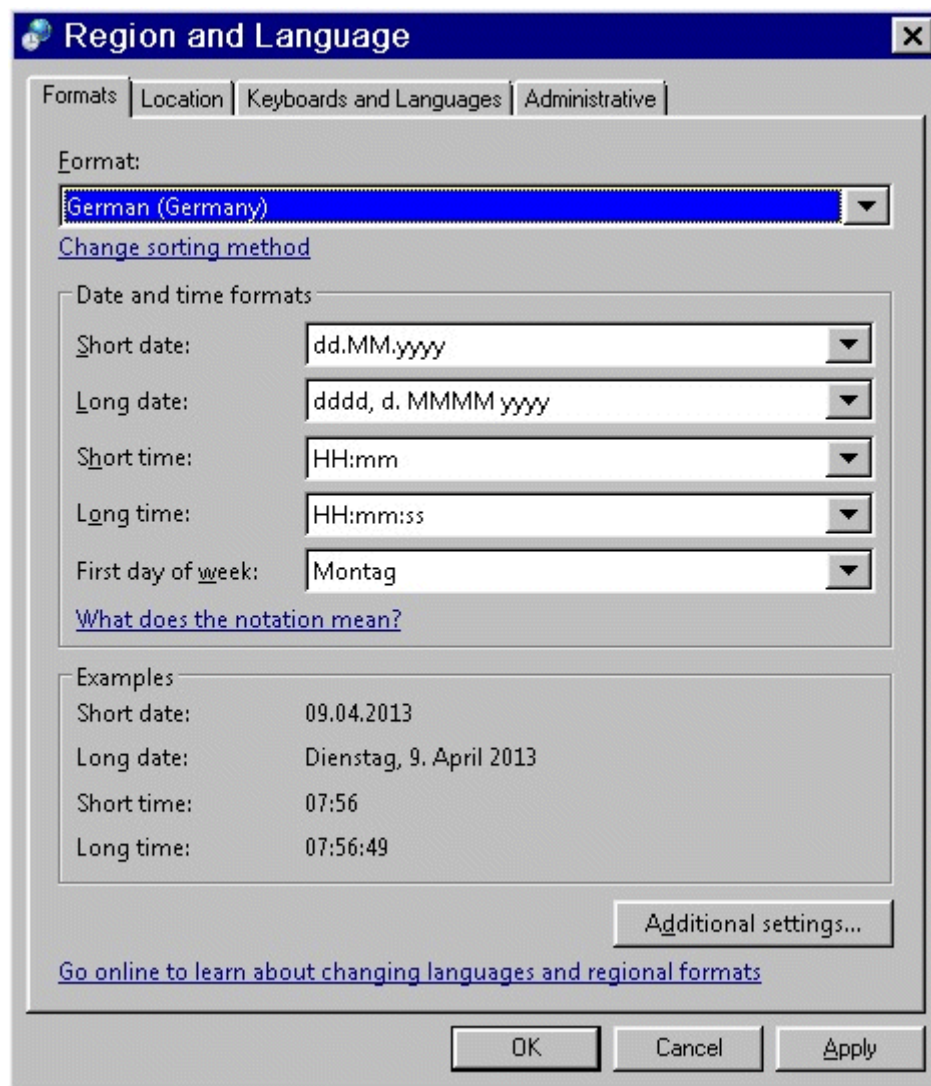

This rule is now with VE11 also valid for Japanese and Chinese language.

4.1.8 Language settings

Supported languages for meduser


- English (US)
- German (Germany)
- French (France)
- Spanish (Spain)
- Japanese

- Chinese (PRC)

4.1.8.1 Keyboard and Language Settings - Japanese


- In the Service UI, select Options -> Configuration
- Select the DICOM Character Set > in the drop-down menu DICOM default character set, select **Japanese, Katakana**

Fig. 100: Configuration - DICOM default character set

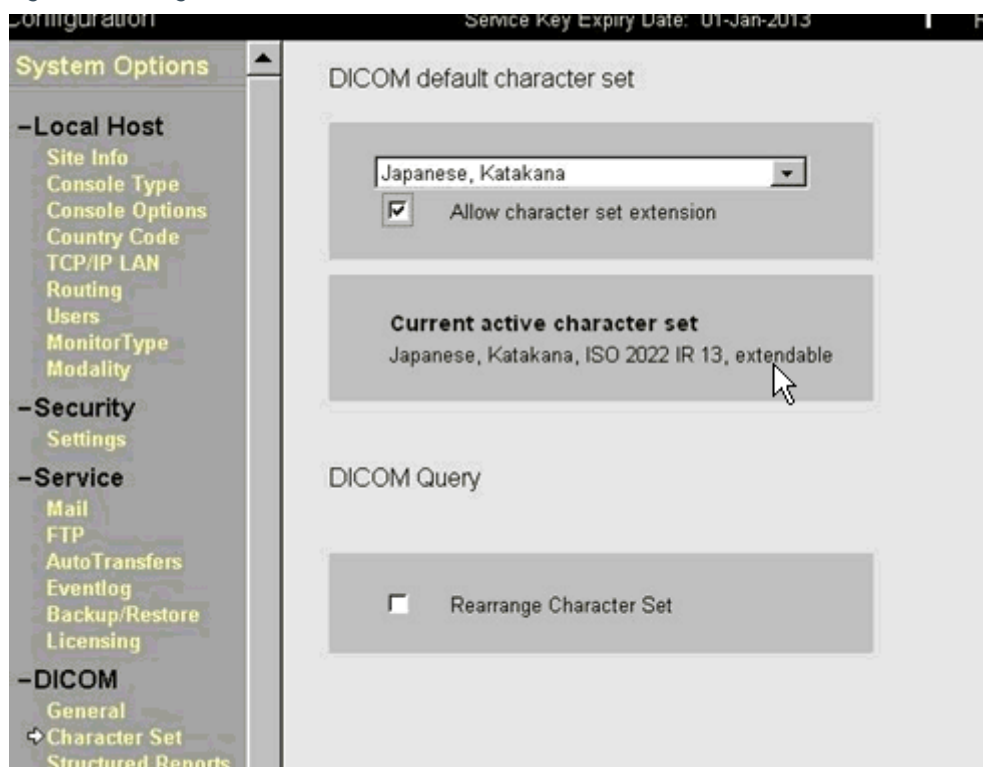


- Click **OK** to confirm the change in character set

Fig. 101: Configuration - character set data successfully saved

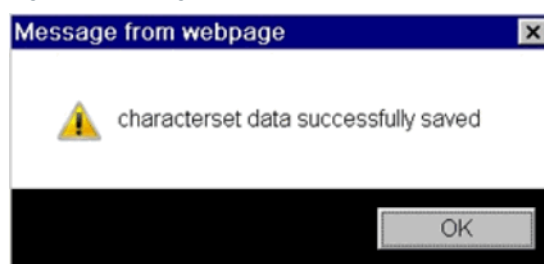


- Click **OK** to restart the application

Fig. 102: Configuration - a restart of the application SW is necessary to make changes valid

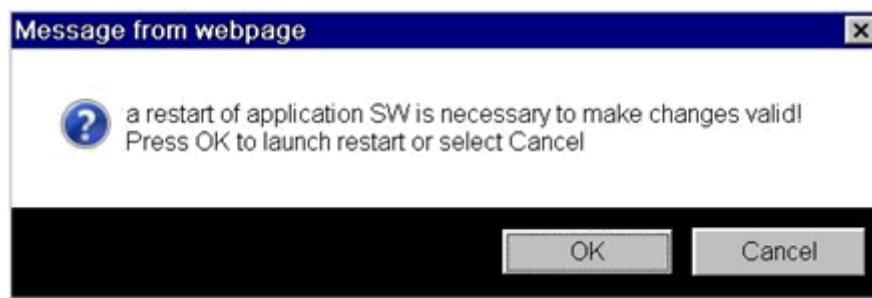


- Select **Options > Configuration**
- Select **Regional and Language Options**.

Fig. 103: Configuration panel

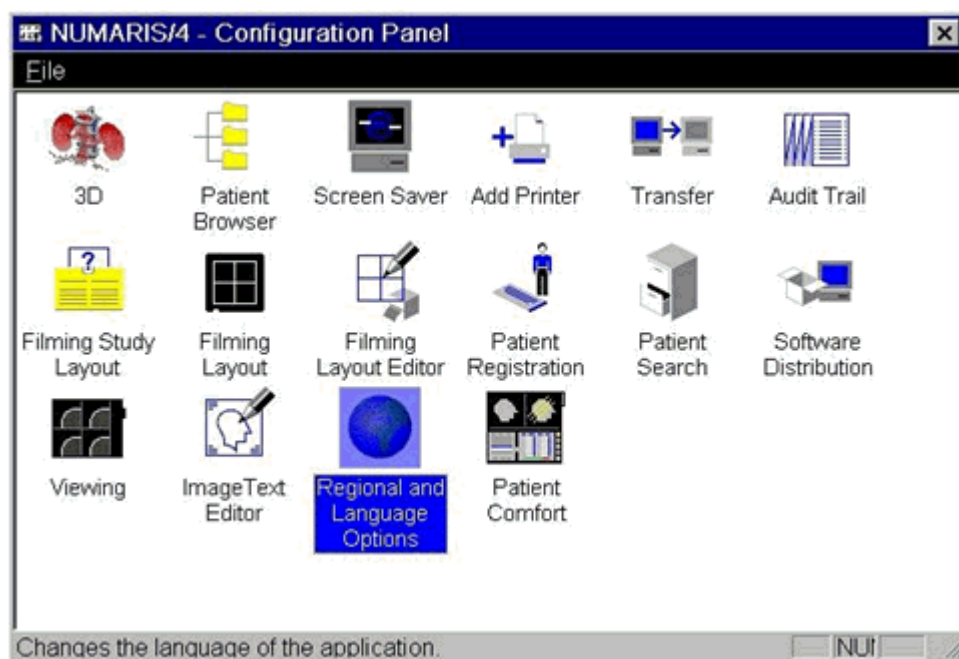

- Select **Japanese (Japan)** in the drop-down menu below "Format"
- Click **Additional settings**

Fig. 104: Regional and Language jp

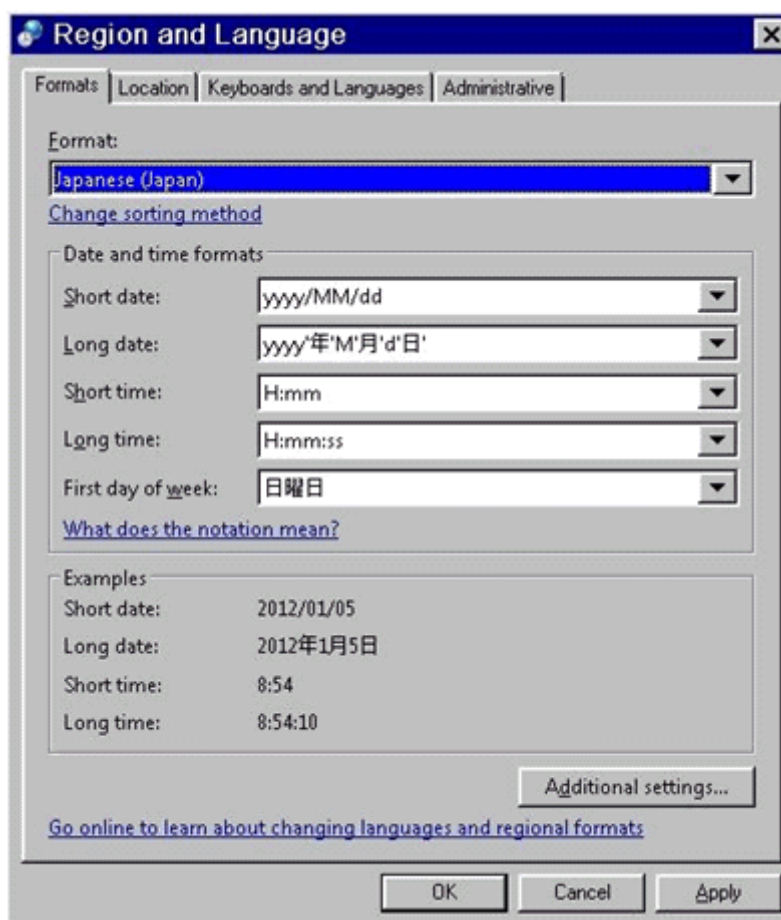


- Select the **Date** tab card
- Select the Calendar drop-down menu

Fig. 105: Calendar type Japan


和暦

Fig. 106: Customize Format - Calendar type

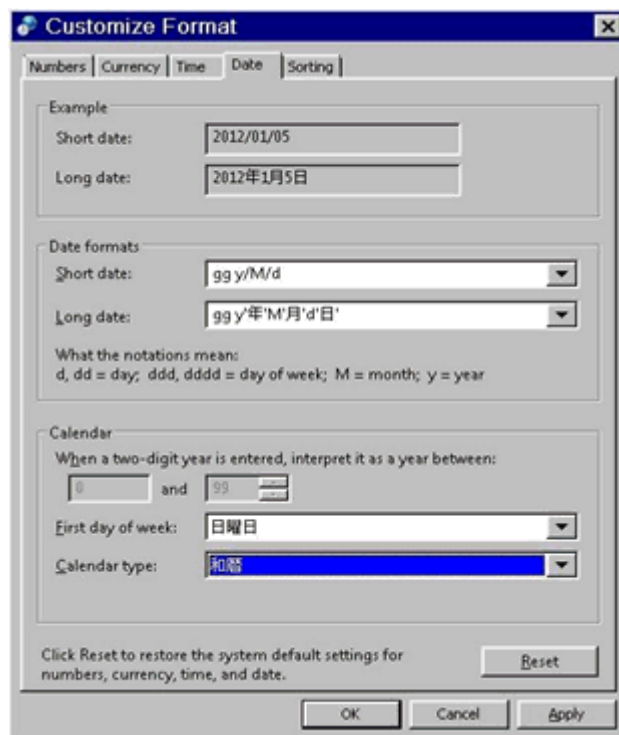

- In the drop-down menu "Current location:" select **Japan**
- Click **OK**

Fig. 107: Region and Language - Location Japan

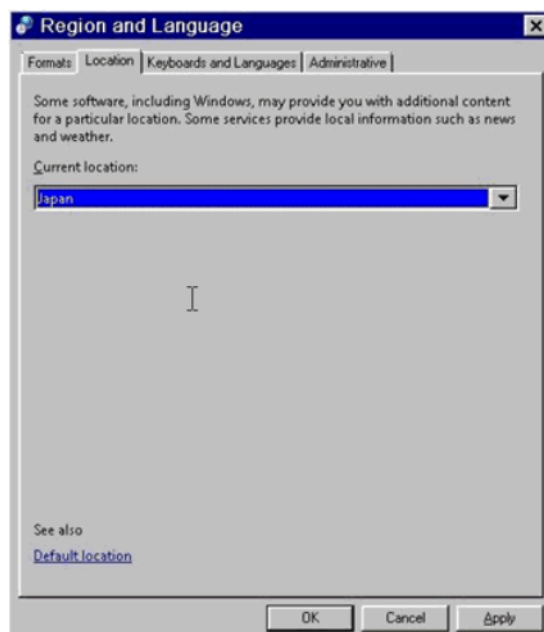


- Select the **Keyboards and Languages** tab card
- In the drop-down menu "Choose a display language:" select:

Fig. 108: Display language Japanese


日本語

Fig. 109: Region and Language - Keyboards and Languages Japan

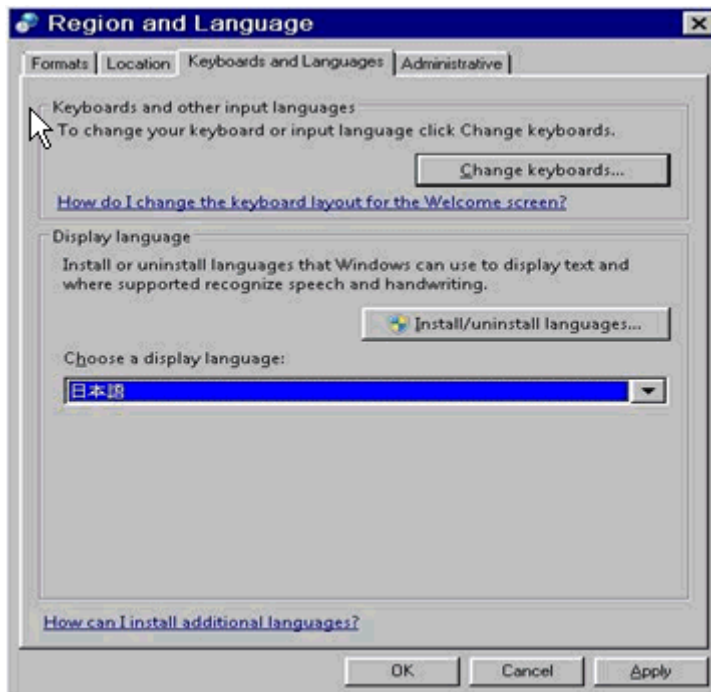


- Select the **General** tab card
- In the drop-down menu "Default input language:" select **Japanese (Japan) - Microsoft IME**:
- Add the Japanese keyboard if not available
- Select **Japanese (Japan) Keyboard Microsoft IME**

Fig. 110: Text Services and Input Languages - General Japan

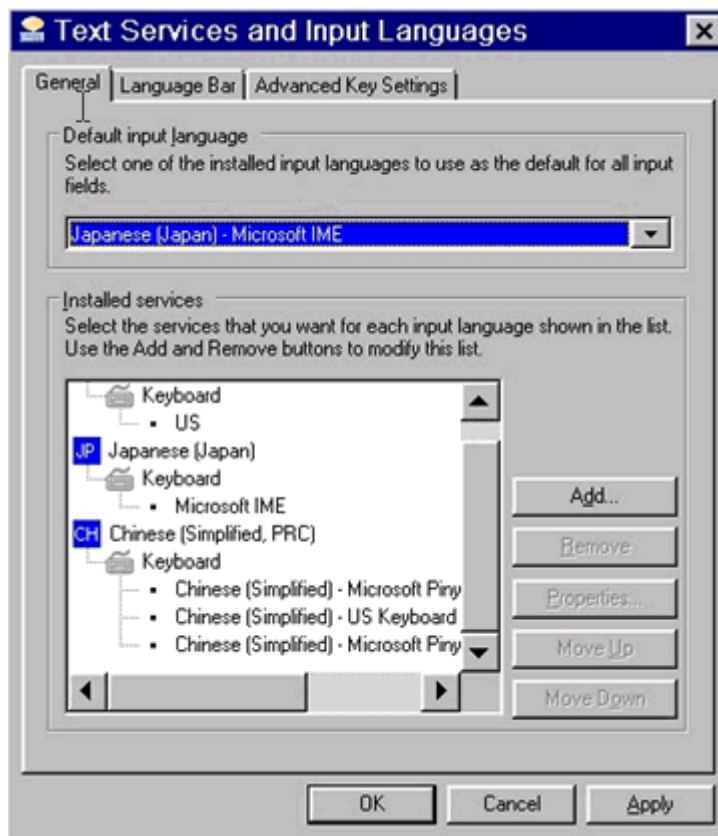

- Set the checkbox **IME** and click **OK**

Fig. 111: Add Input Language Japanese

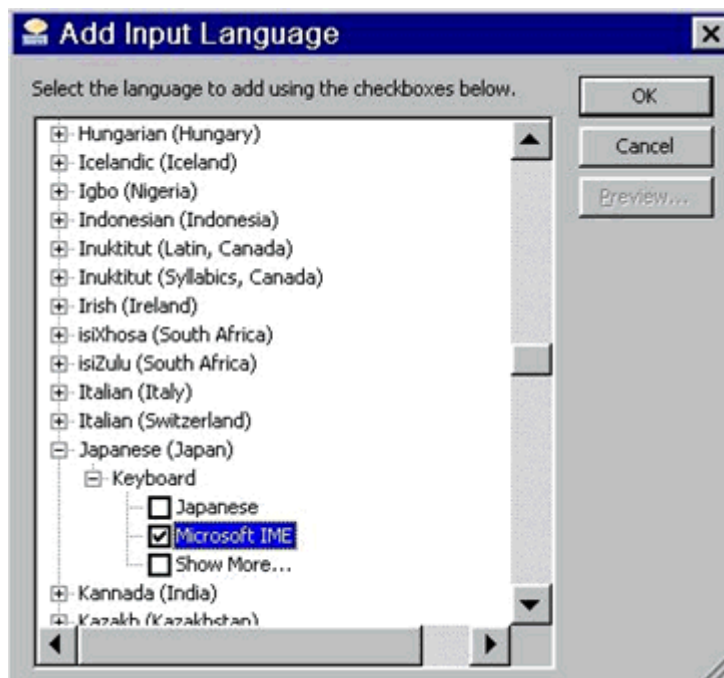

- Select the **Language Bar** tab card
- Select the **Floating On Desktop** radio button


- Set the checkboxes **Show additional Language bar icons in the taskbar** and **Show text labels on the Language bar**
- Click **OK**

Fig. 112: Text Services and Input Languages



- Select the **Administrative** tab card
- Click **Change system locale**

Fig. 113: Region and Languages - Administrative

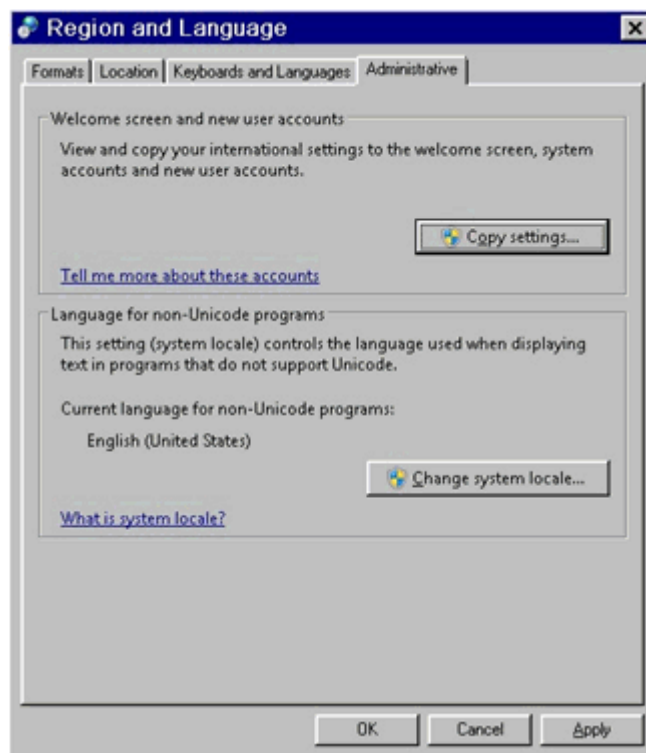

- Click **Apply** in the popup window

Fig. 114: Change Regional Options




- Select Japanese (Japan) in the next popup window
- Click **OK**

Fig. 115: Region and Language Settings - Japanese




- In the following popup window select **Restart now**

Fig. 116: Change System Locale - Restart now

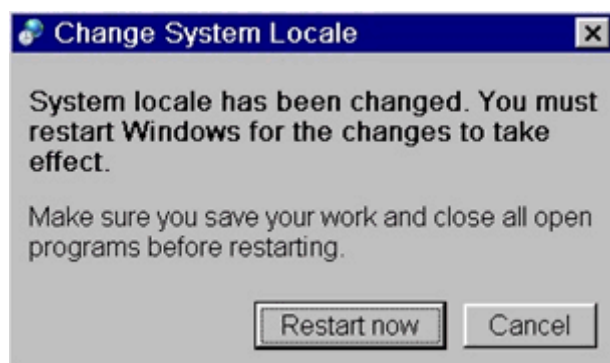


- Confirm this popup with **OK**

Fig. 117: OK - jp

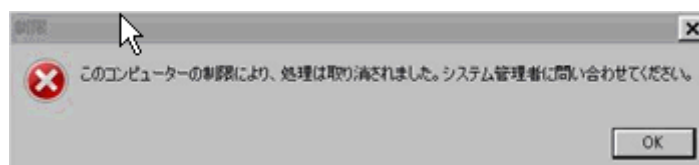

- Confirm the warning with **OK**


Fig. 118: Warning OK - jp

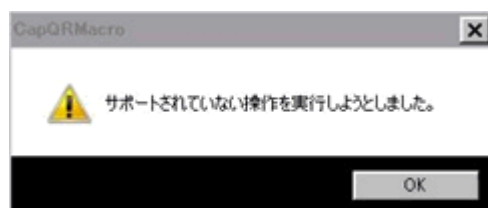

4.1.8.2 Switch to English from Japanese


- Select **Options > Service > Local Service**

Fig. 119: Options - Service - Local Service.jp

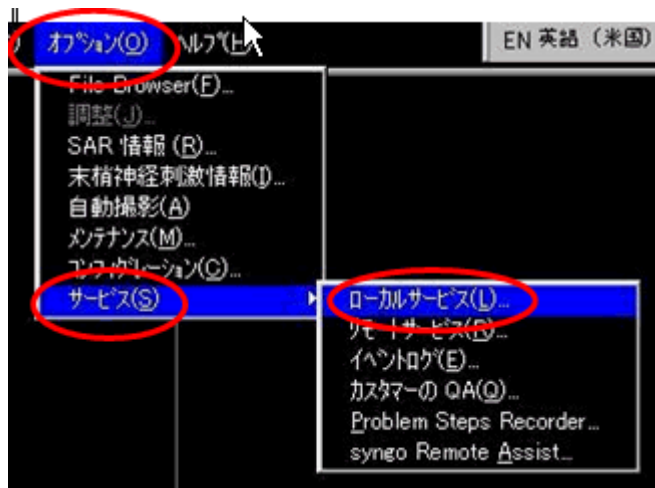

- Select **DICOM Character Set >** in the drop-down menu "DICOM default character set" select **Latin Alphabet No. 1**

Fig. 120: Configuration - Character Set - DICOM default character set (Latin Alphabet No 1) .jp

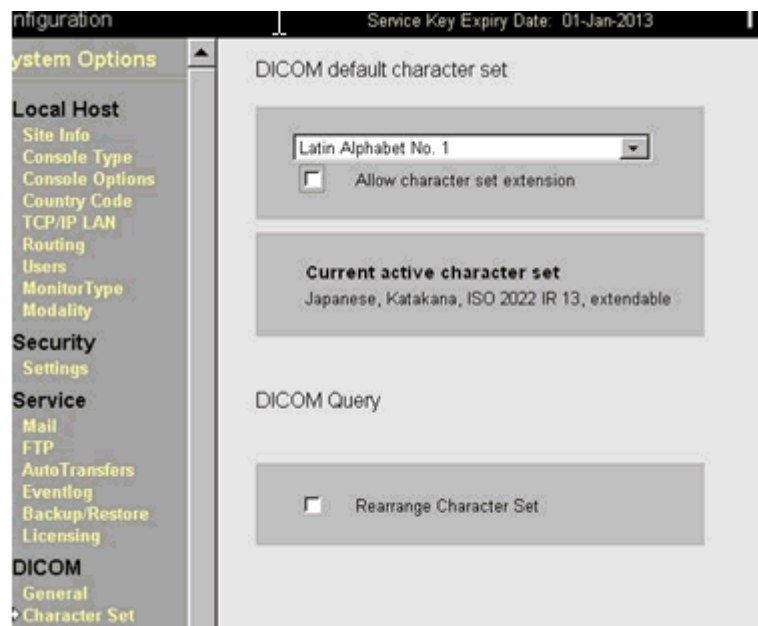


- Confirm "character set data successfully saved" popup with **OK**

Fig. 121: Character set data successfully saved - OK



- Confirm the "application restart popup" with **OK**

Fig. 122: A restart of the application SW is necessary - OK.jp

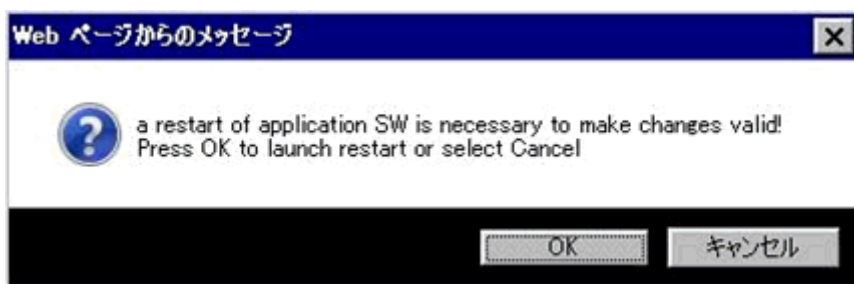


- When the application is up and running again, close the windows "Personalization" and "Local Service"

Fig. 123: Personalization and Local Service - jp




- Click Options > Configuration

Fig. 124: Options Configuration - jp

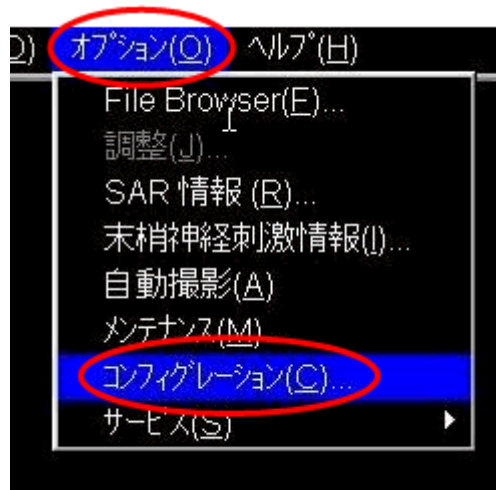

- In the Numaris/4 configuration panel, double-click **Regional and Language Options**


Fig. 125: Regional and Language Options - jp

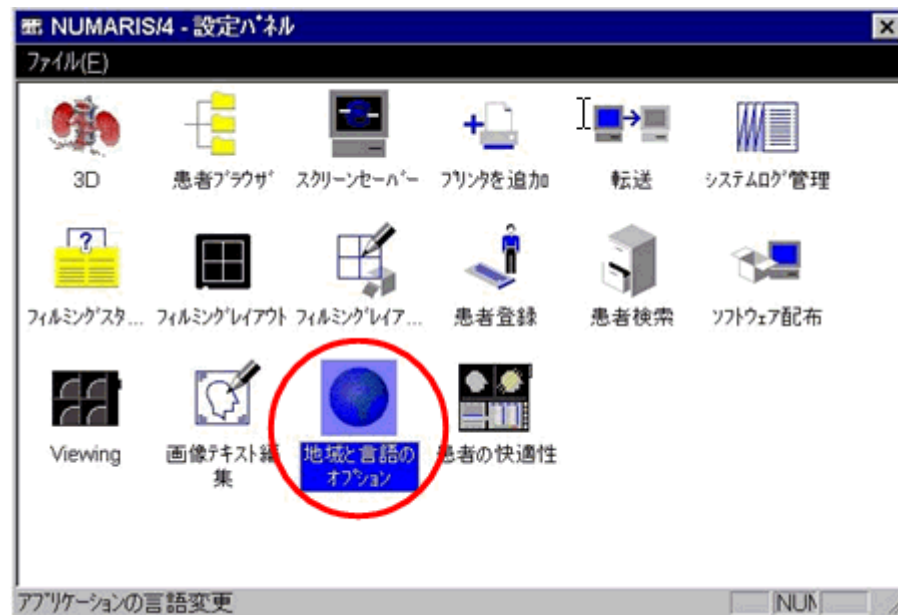

- Set the date time format as follows

Fig. 126: Region and Language - jp


英語 (米国)

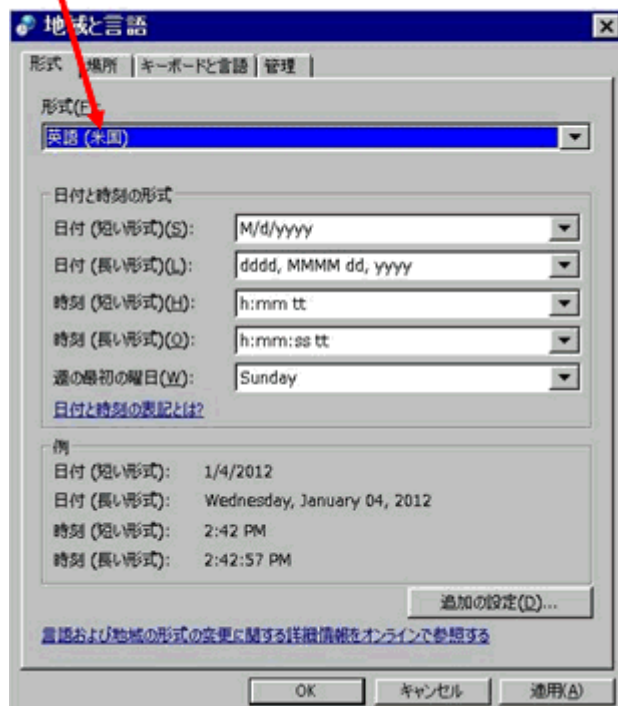

- Select the second tab
- Click OK


Fig. 127: Region and Language - Location jp

米国

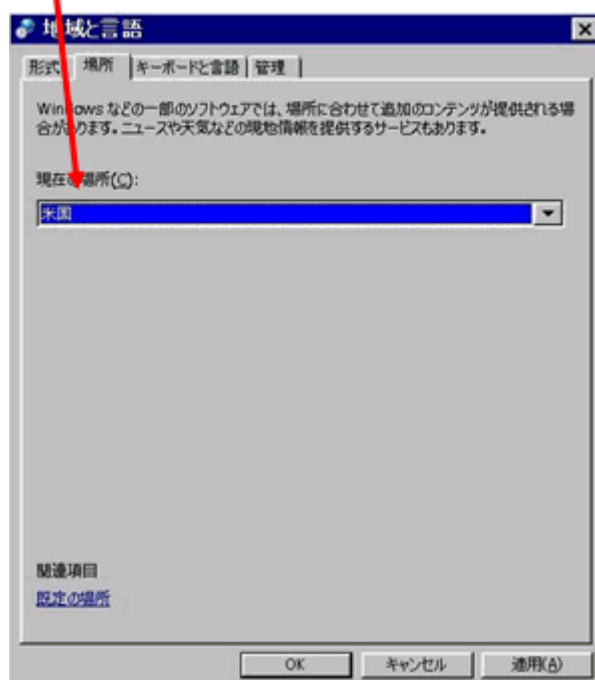

- Select English


Fig. 128: Region and Language - Keyboards and Languages jp

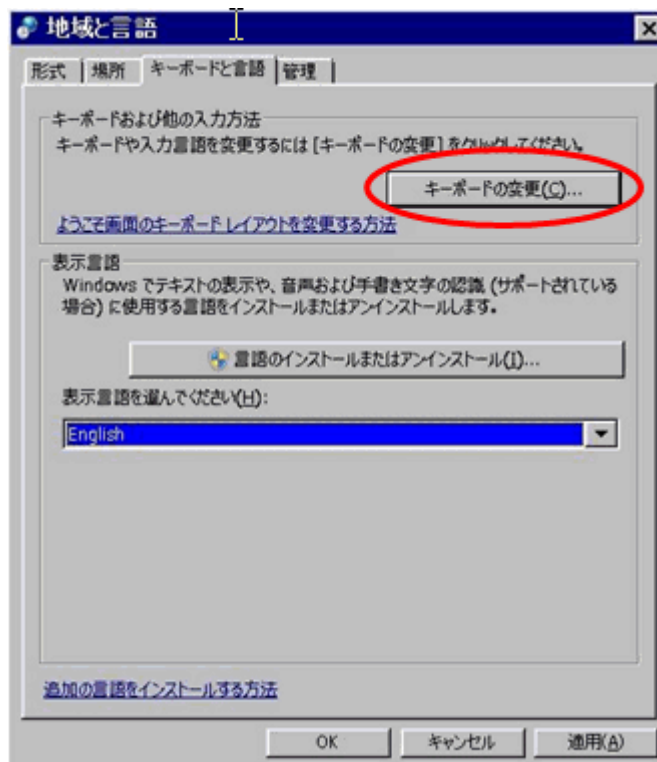

- Click the first tab from the left
- Select the following in the left drop-down menu

Fig. 129: Text Services and Input Languages jp


英語 (米国) - US

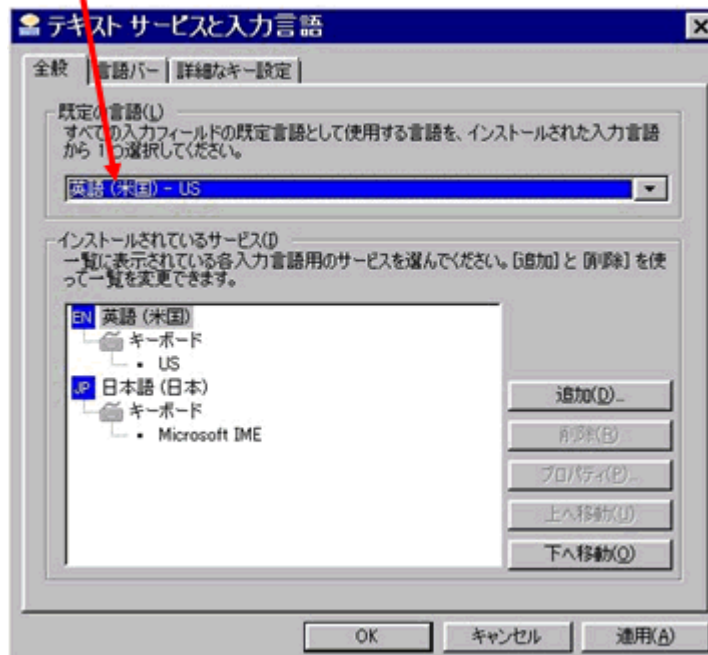

- Select the second tab > select the mid radio button

- Click OK

Fig. 130: Text Services and Input Languages jp

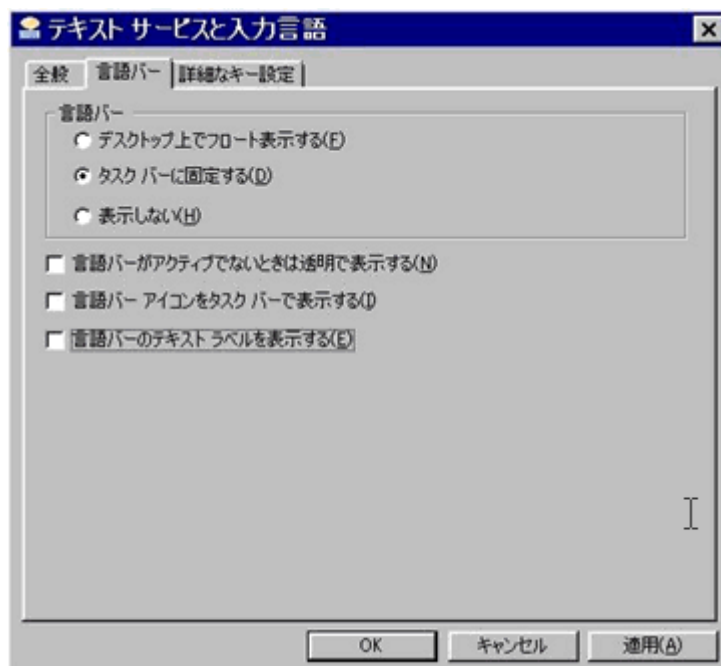

- Select the fourth tab > select the button

Fig. 131: Region and Language - Change system locale jp (zoomed)



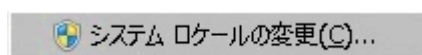


Fig. 132: Region and Language - Change system locale.jp

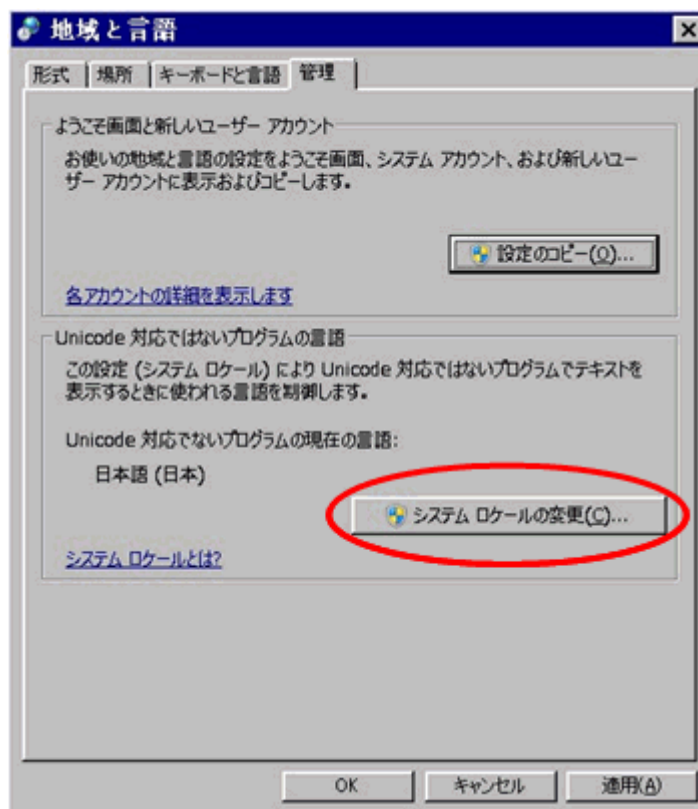


- Click **Apply**.

Fig. 133: Change Regional Options - Apply.jp

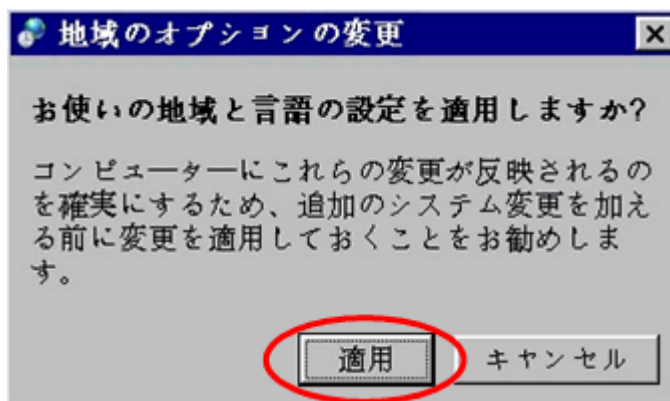


- Select the drop-down menu as follows and click **OK**

Fig. 134: Region and Language Settings - US English jp

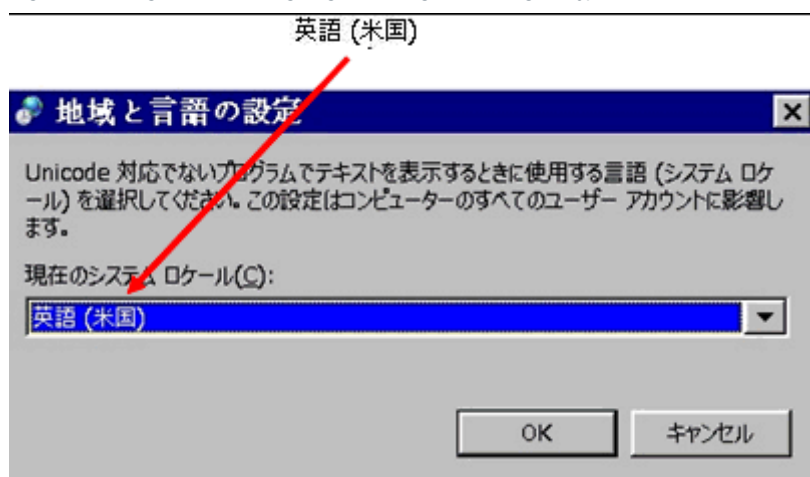

- Proceed as follows:

Fig. 135: Change System Locale - Restart now ch

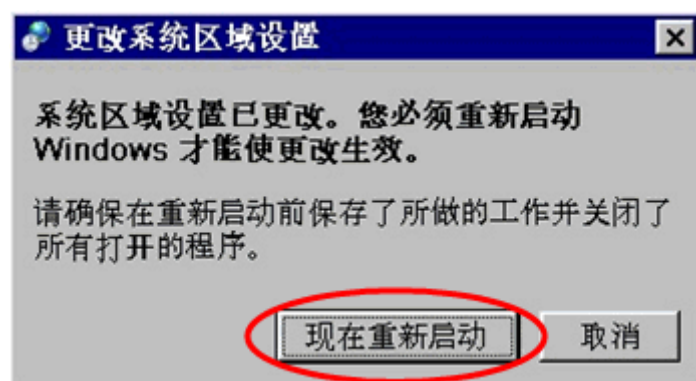

- Confirm the restriction popup with OK


Fig. 136: Confirm Restriction popup with OK



4.1.8.3 Keyboard and Language Settings - Chinese


- In the Service UI, select **Options** -> **Configuration**
- Select the DICOM **Character Set** > in the drop-down menu "DICOM default character set" select **Chinese**

Fig. 137: Configuration - DICOM default character set

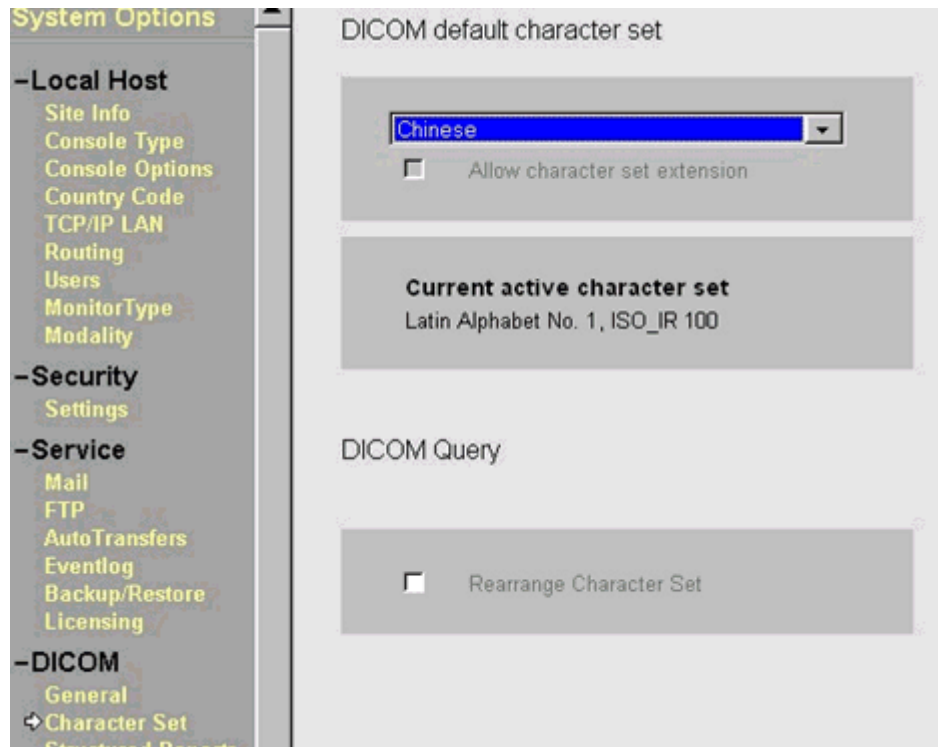


- Click **OK** to confirm the change in character set

Fig. 138: Configuration - character set data successfully saved

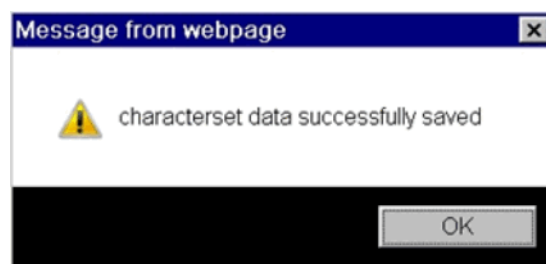


- Click **OK** to restart the application

Fig. 139: Configuration - a restart of the application SW is necessary to make changes valid

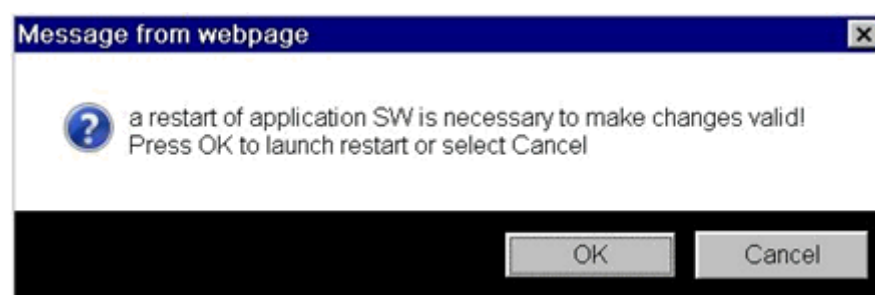


- Select **Options > Configuration**
- Select **Regional and Language Options**.

Fig. 140: Configuration panel

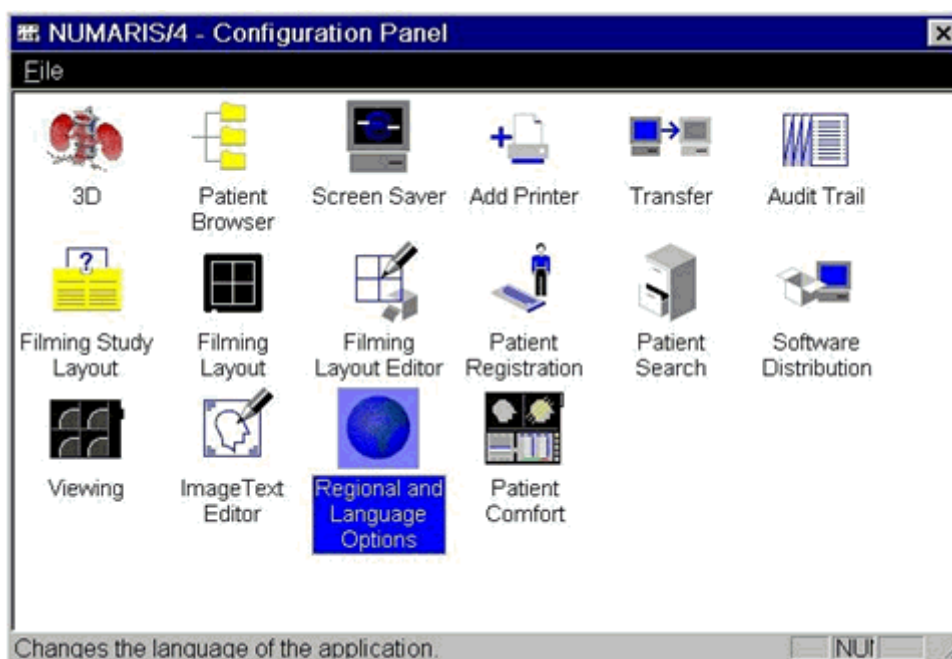


- Select **Chinese (Simplified, PRC)** in the drop-down menu below "Format"

Fig. 141: Regional and Language

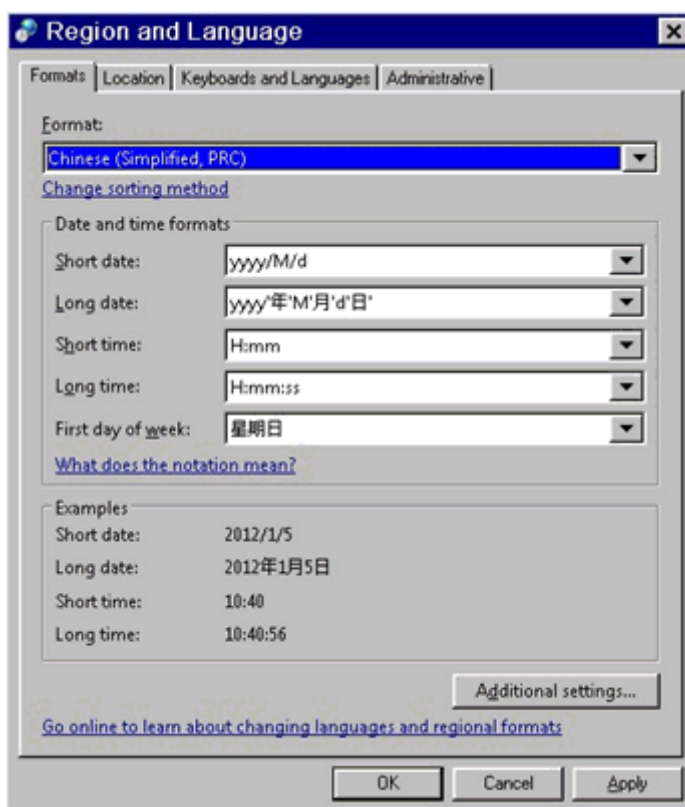

- Select the **Location** tab
- Select **Chinese** in the drop-down menu

Fig. 142: Region and Language - Location China

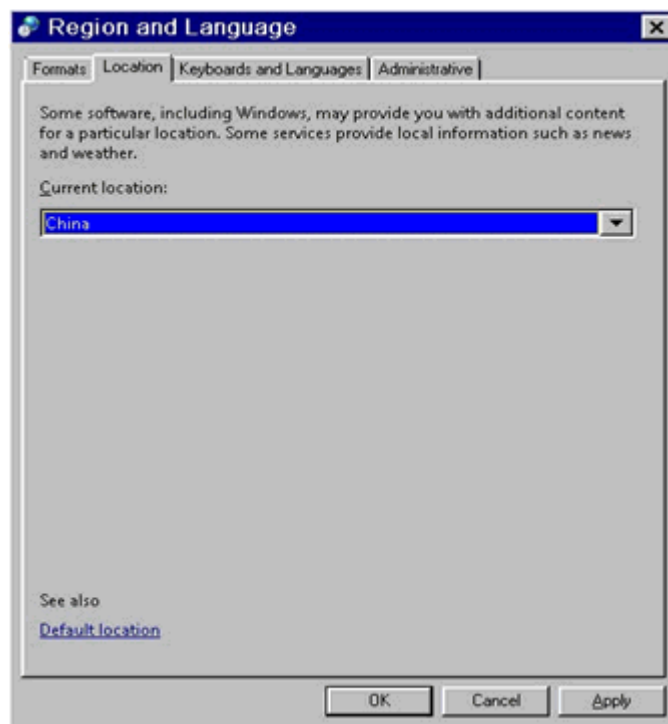


- Select the **Keyboards and Languages** tab
- Select

Fig. 143: Choose display language - Chinese

中文(简体)


- in the drop-down menu **Choose a display language** tab

Fig. 144: Region and Language - Keyboards and Languages

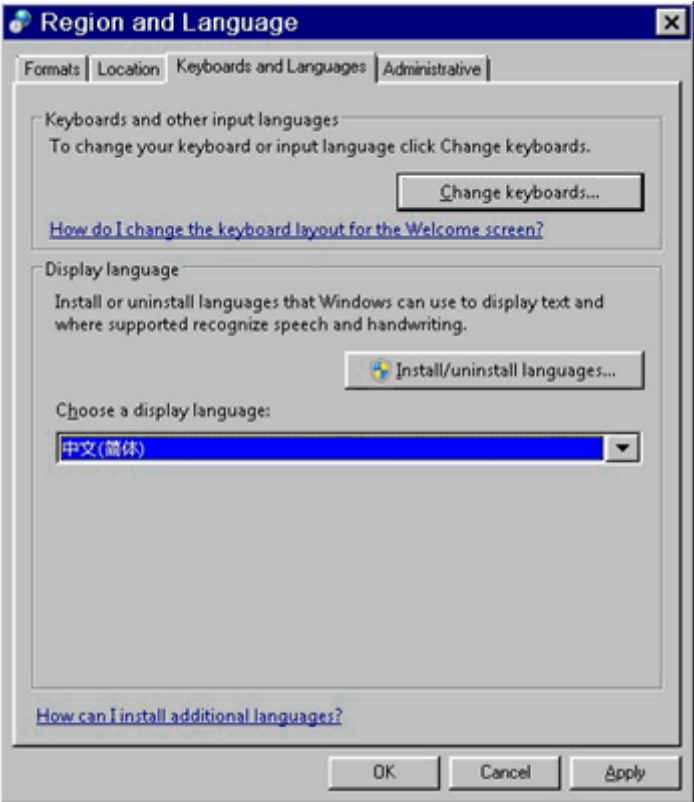


- Select the **General** tab.
- Select the drop-down menu **Chinese (Simplified, PRC) - Chinese (Simplified) - Microsoft Pinyin ABC Input Style**

Fig. 145: Text Services and Input Languages - General

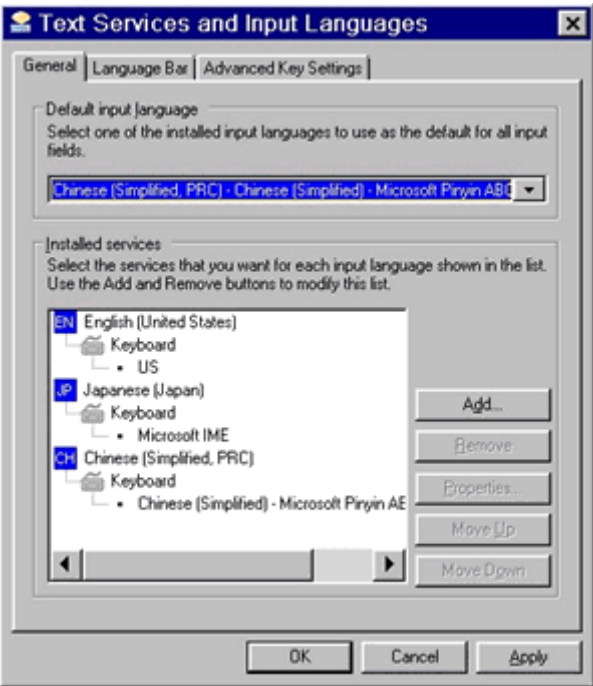

- Set the checkbox “Chinese (Simplified, PRC) - Chinese (Simplified) - Microsoft Pinyin ABC Input Style”

Fig. 146: Add Input Language

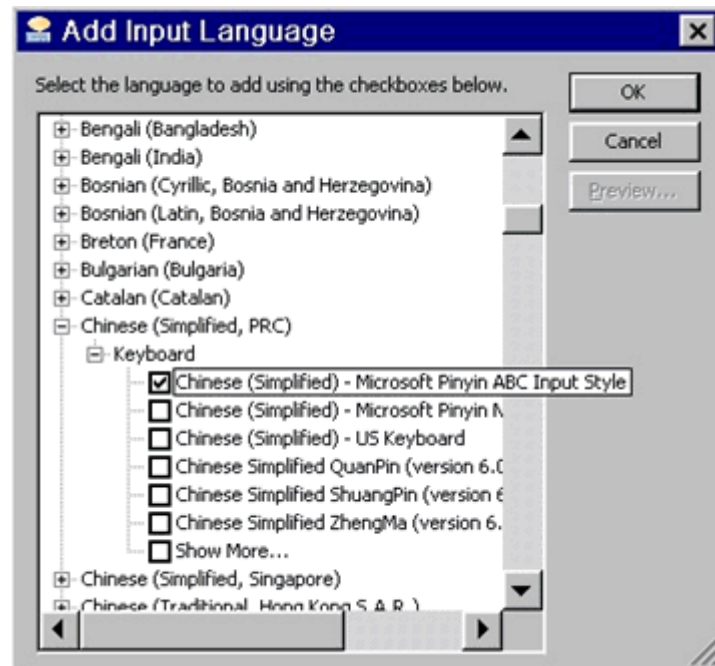


- Select the **Language Bar** tab card
- Select the **Floating On Desktop** radio button
- Set the checkbox **Show additional language bar icons in the task bar**
- Set the checkbox **Show text labels on the language bar**

Fig. 147: Text Services and Input Languages




- Select the **Administrative** tab
- Click **Change system locale**

Fig. 148: Region and Languages - Administrative

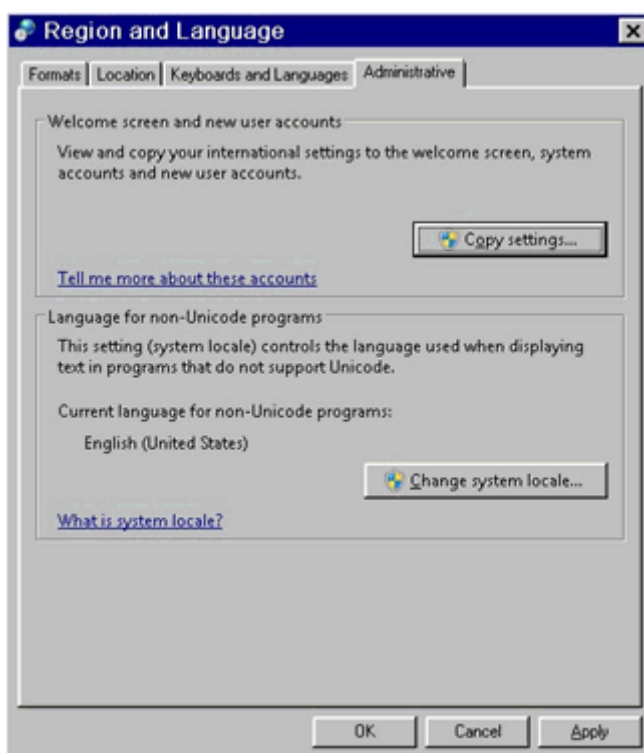


- Click **Apply** to confirm the region and language changes

Fig. 149: Change Regional Options

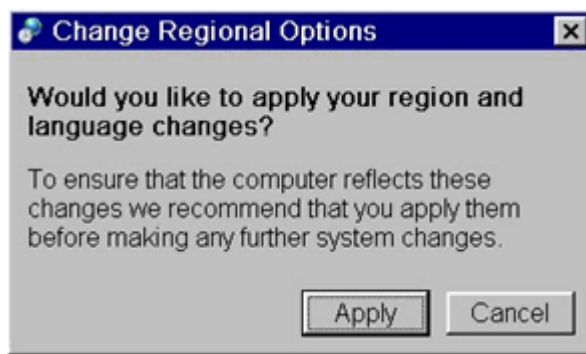


- Select **Chinese (Simplified, PRC)** in the drop-down menu "Current system locale"

Fig. 150: Region and Language Settings - Chinese

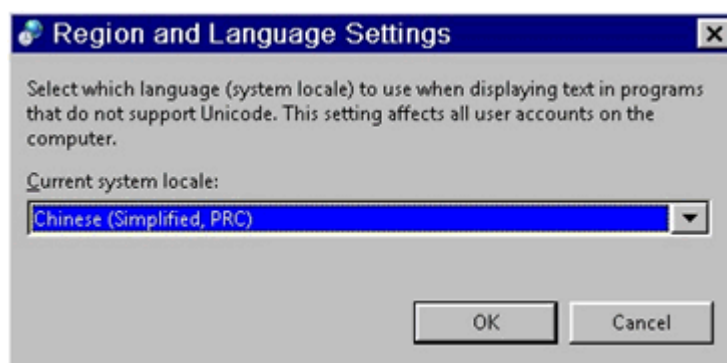


- Click **Restart now**

Fig. 151: Change System Locale - Restart now

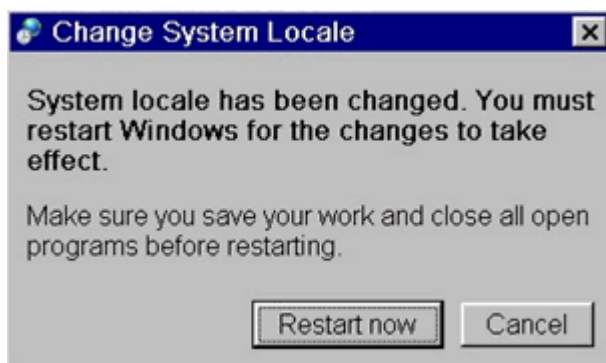

- Log in as meduser, wait until the application is up and running
- Confirm the following popups with **OK**

Fig. 152: OK - ch

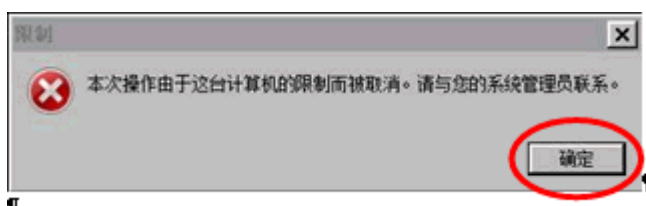

- Confirm the warning with **OK**


Fig. 153: Warning OK - ch



4.1.8.4 Switch to English from Chinese


- In the Service UI, select **Options** -> **Configuration**
- Select the DICOM **Character Set** > in the drop-down menu "DICOM default character set" select **Chinese**
- Open the Local Service service GUI

Fig. 154: Options - Service - Local Service ch

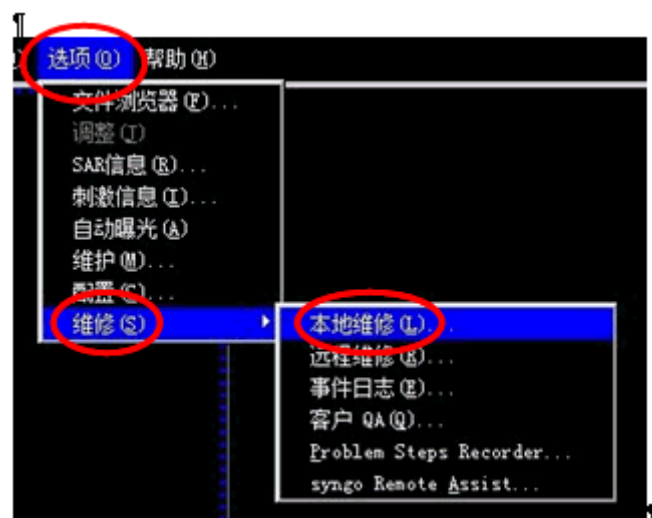


- Set the DICOM Character Set to **Latin Alphabet No. 1**

Fig. 155: Configuration - Character Set - DICOM default character set (Latin Alphabet No 1)

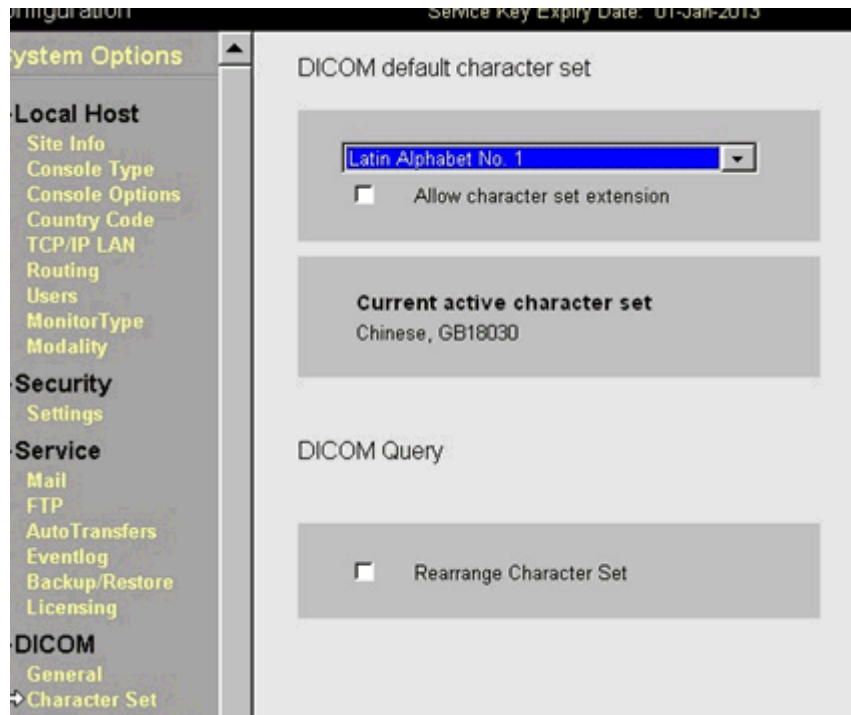


- Confirm the "character set data successfully saved" popup with **OK**

Fig. 156: Character set data successfully saved - OK ch

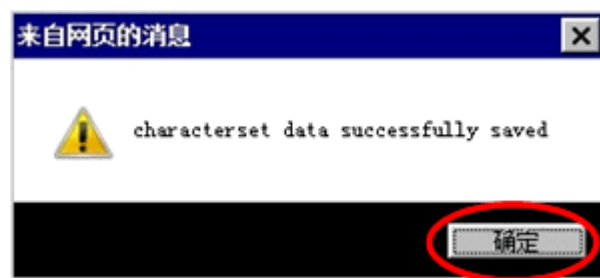

- Confirm the application restart popup with **OK**

Fig. 157: A restart of the application SW is necessary - OK ch

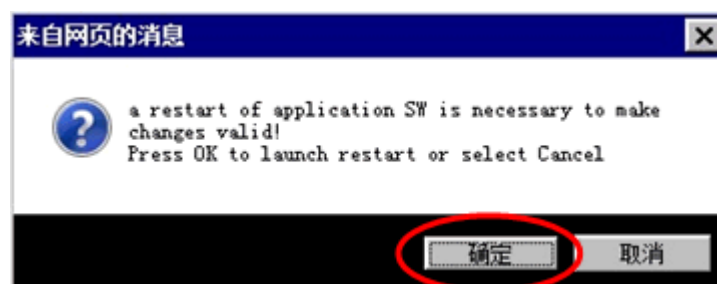


- When the application is up and running again, close the windows "Personalization" and "Local Service"

Fig. 158: Personalization and Local Service - ch



- Click **Options > Configuration**

Fig. 159: Options Configuration - ch

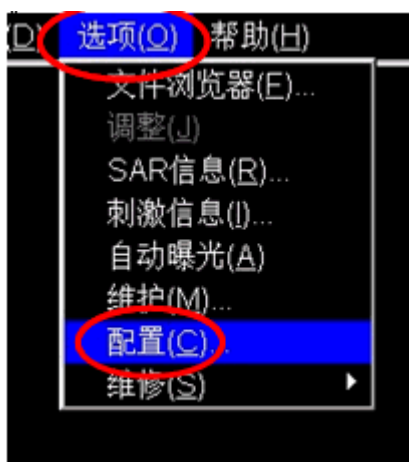


- In the Numaris/4 configuration panel, double-click **Regional and Language Options**

Fig. 160: Regional and Language Options - ch




- Set the date time format as follows

Fig. 161: Region and Language - ch

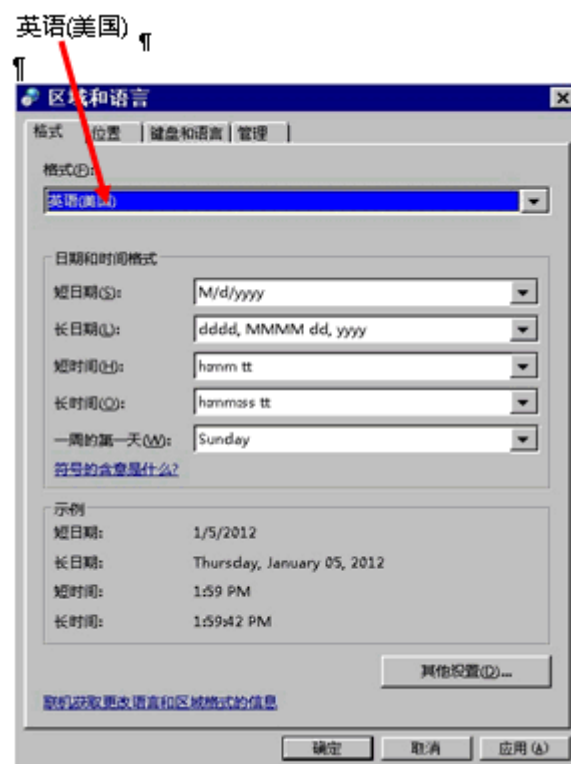


- Select the second tab
- Click OK.

Fig. 162: Region and Language - Location ch

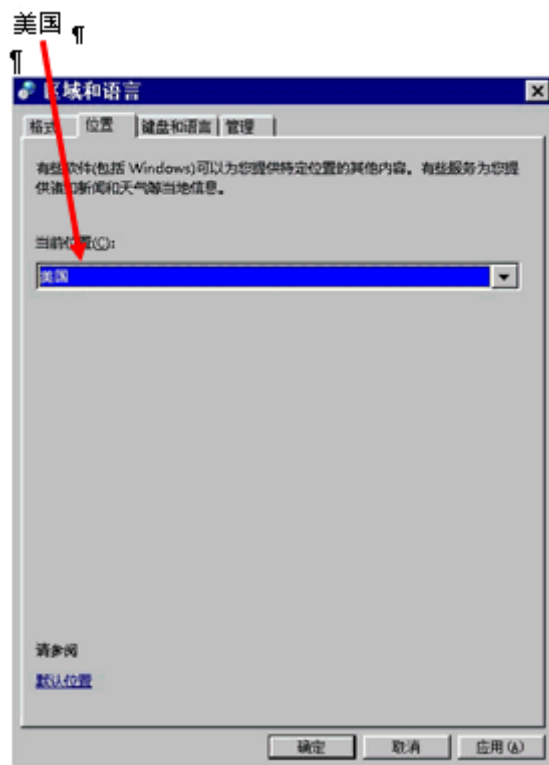


- Select **English**

Fig. 163: Region and Language - Keyboards and Languages ch

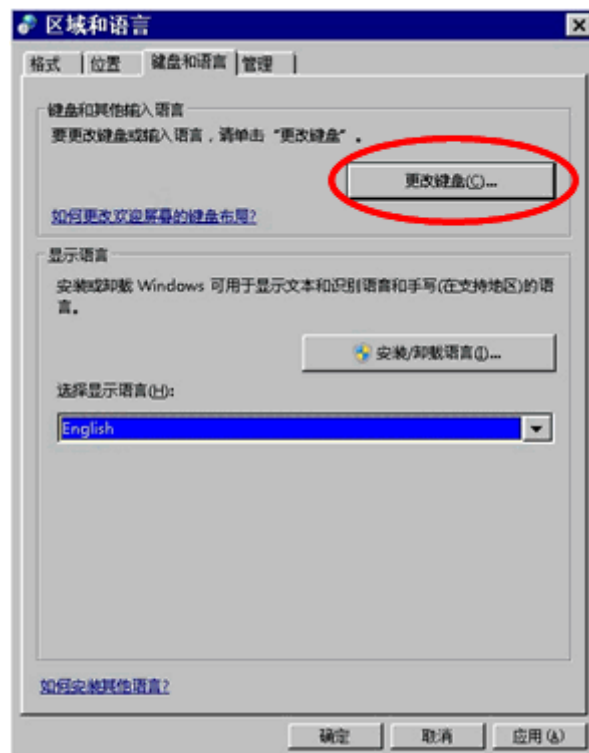


- Click the first tab from the left
- Select the following in the left drop-down menu

Fig. 164: Text Services and Input Languages ch

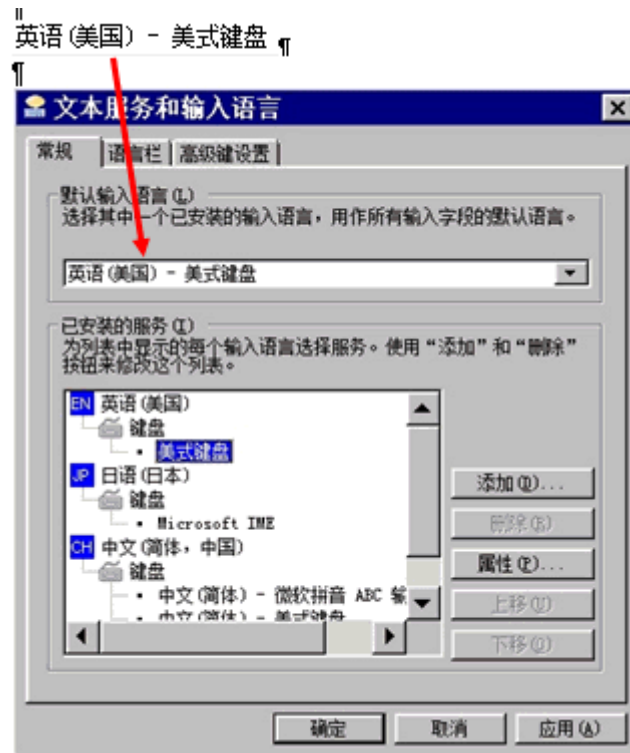

- Select the second tab > select the mid radio button
- Click OK

Fig. 165: Text Services and Input Languages ch

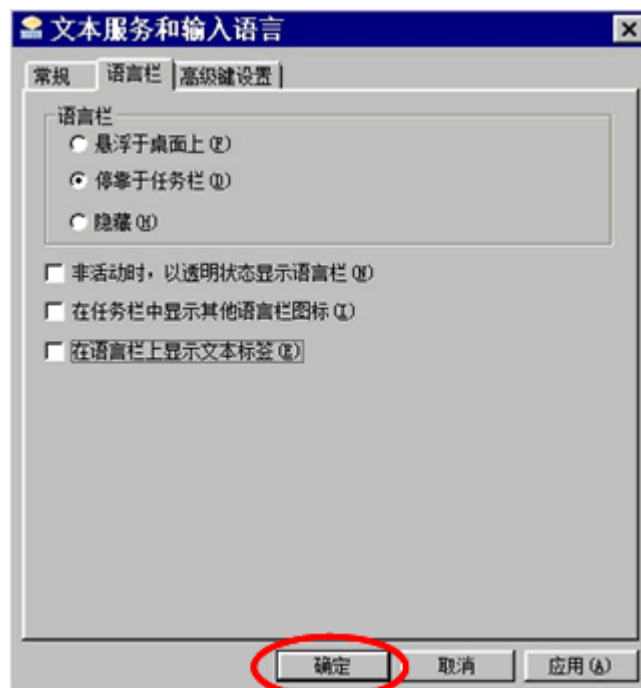

- Select the fourth tab > select the button

Fig. 166: Region and Language - Change system locale ch (zoomed)



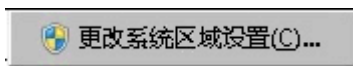


Fig. 167: Region and Language - Change system locale ch

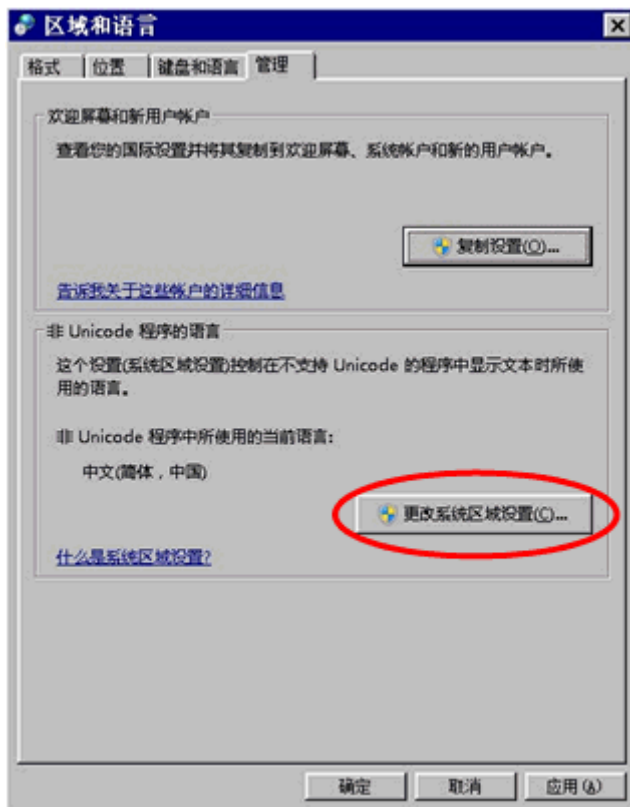


- Click **Apply**.

Fig. 168: Change Regional Options - Apply ch

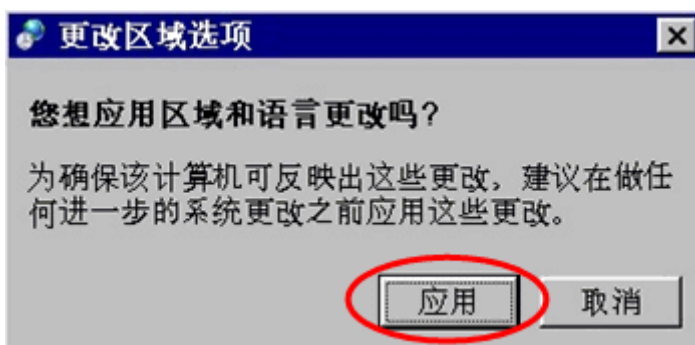


- Select the drop-down menu as follows and click **OK**

Fig. 169: Region and Language Settings - US English ch

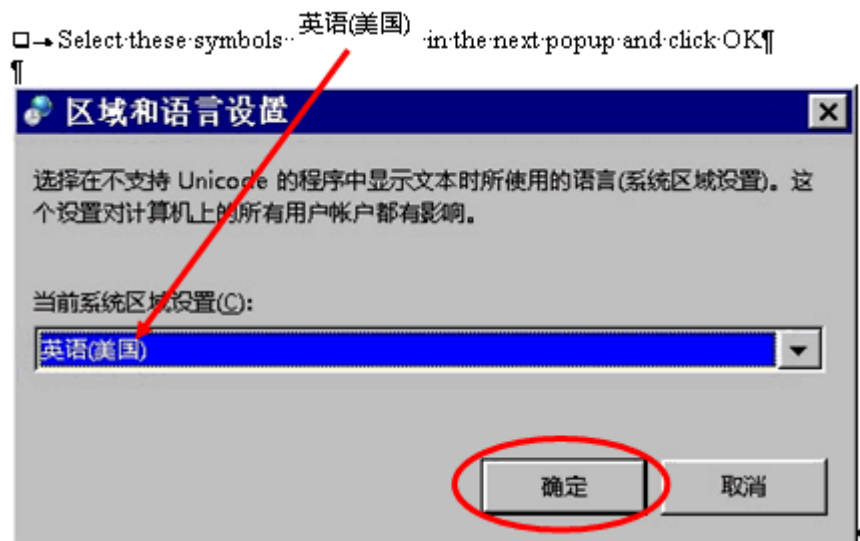

- Proceed as follows:

Fig. 170: Change System Locale - Restart now ch

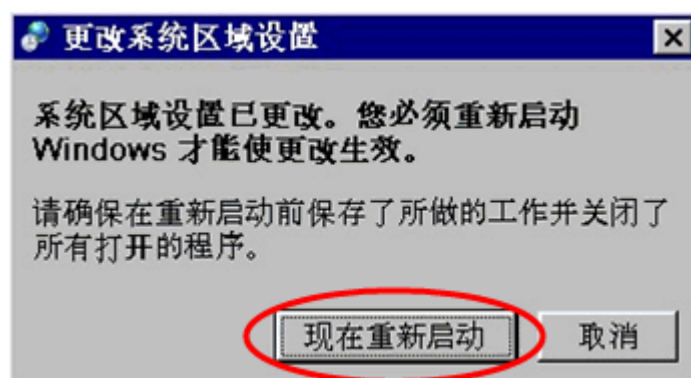

- Confirm the restriction popup with OK

Fig. 171: Confirm Restriction popup with OK



4.1.9 System management installation workaround

4.1.9.1 Description:

If system management under Agent Controls is deinstalled, if the checkboxes

- SD -Asset - Management Agent
- Event Management Agent

are not marked (→ Fig. 172 Page 140), or if the installation is started again (both the checkboxes are marked), the Event Management Agent installation fails (→ Fig. 173 Page 141).

Fig. 172: MNP deinstalled


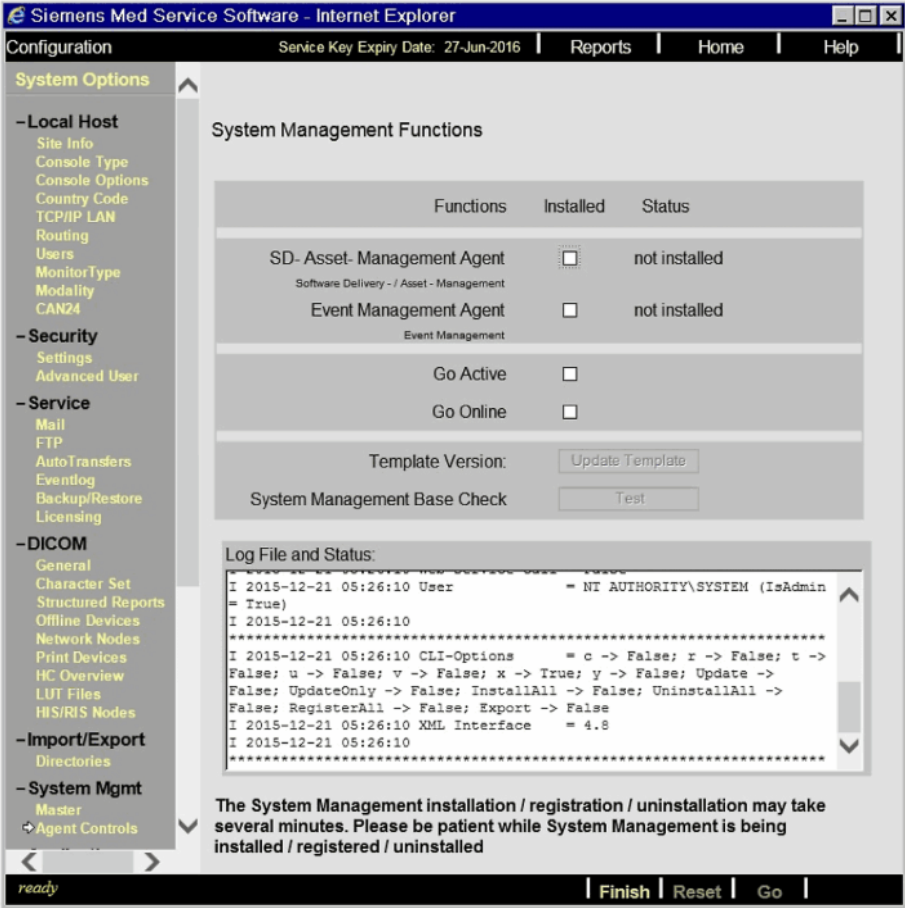


Fig. 173: MNP Event Agent installation failed



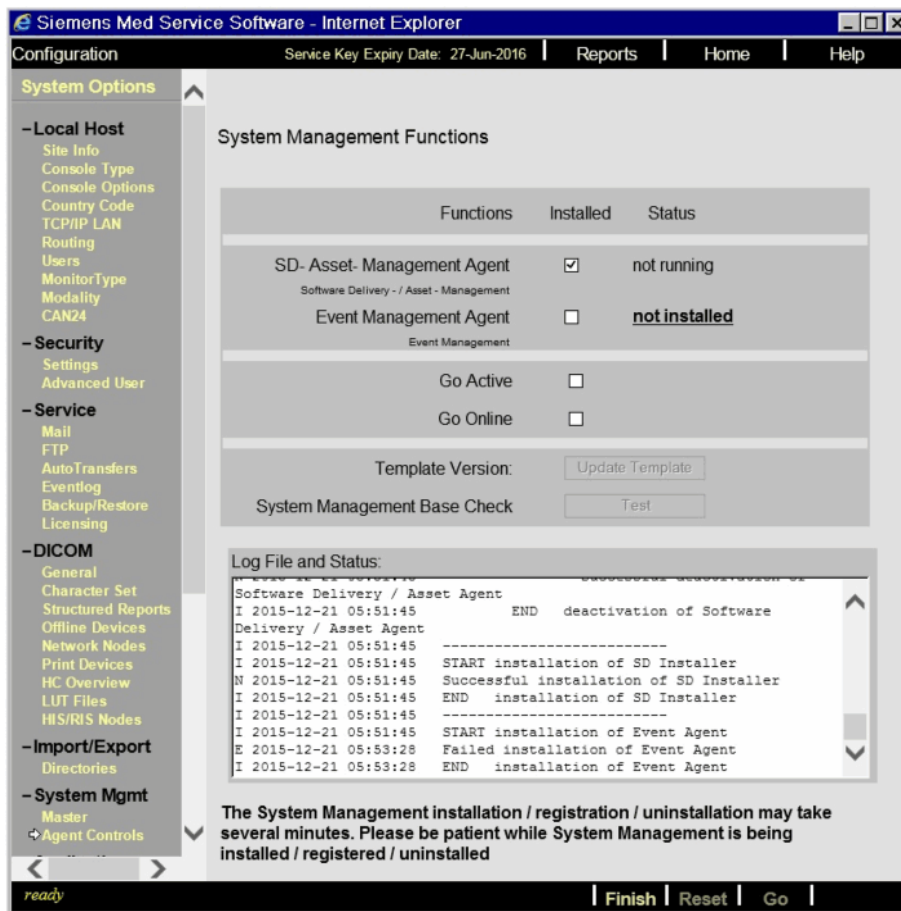



4.1.9.2 Solution:


The following procedure describes the workaround for this.

- Set Embedded Control in update mode

ACTION	MRAWP/MRWP
■ Close the local service UI.	☐
■ Enable Advanced User	☐
■ Start the MREmbeddedControlGUI from the Start menu Start ⇒ All Programs ⇒ MREmbeddedControl ⇒ MREmbeddedControlGUI	☐

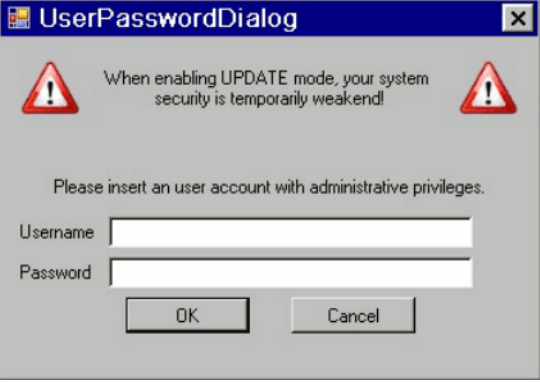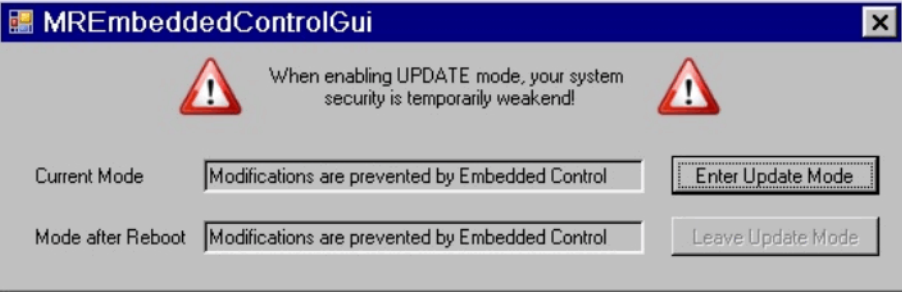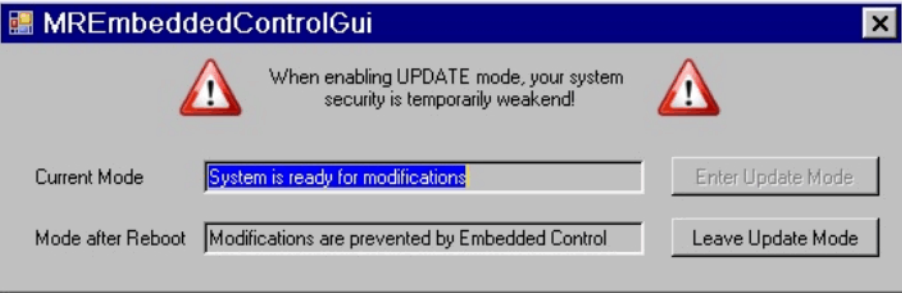ACTION	MRAWP/MRWP
In the popup displayed, the user account with administrative privileges is requested. Fig. 174: EC User Password request  Type in as user "syadmin" and enter the password.	☐
Now the MREmbeddedControlGUI is displayed Fig. 175: MREmbeddedControlGUI 	☐
Press the button Enter Update Mode to switch embedded control into the update mode.	☐
Confirm the displayed popup with OK	☐
Now switch embedded control is in update mode. Fig. 176: MREmbeddedControlGUI update mode 	☐
Leave the MREmbeddedControlGUI open, don't close it.	☐
System Management installation and activation	
ACTION	MRAWP/MRWP
Select Options ⇒ Service ⇒ Local Service.	☐

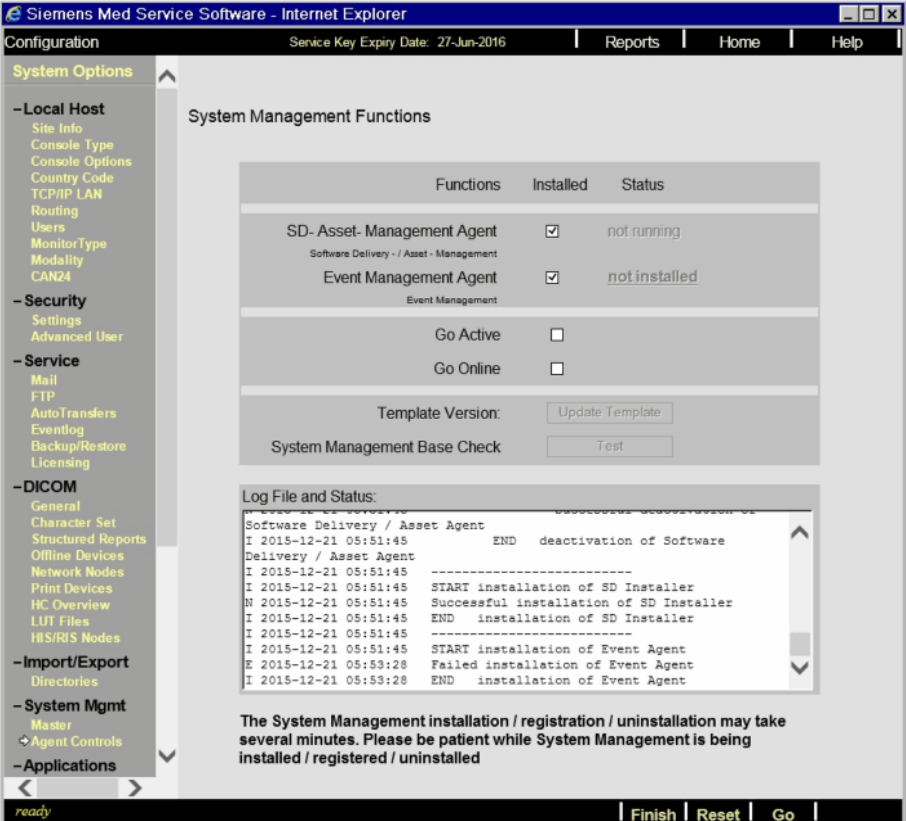ACTION	MRAWP/MRWP															
Enter the valid service key.	☐															
In the local service UI, select Configuration.	☐															
Under System Mgmt. Agent Controls, set the following configuration: Mark the checkboxes "SD- Asset- Management Agent" and "Event Management Agent". Fig. 177: MNP installation new  The screenshot shows the 'Configuration' window of the Siemens Med Service Software. The left sidebar has 'Agent Controls' selected under 'System Mgmt'. The main panel displays 'System Management Functions' with a table: <table><tr><th>Functions</th><th>Installed</th><th>Status</th></tr><tr><td>SD- Asset- Management Agent Software Delivery - / Asset - Management</td><td>☒</td><td>not running</td></tr><tr><td>Event Management Agent Event Management</td><td>☒</td><td>not installed</td></tr><tr><td>Go Active</td><td>☐</td><td></td></tr><tr><td>Go Online</td><td>☐</td><td></td></tr></table> Below the table, there are buttons for 'Update Template' and 'Test'. A 'Log File and Status' section shows the following log entries: <pre>I 2015-12-21 05:51:45 END deactivation of Software I 2015-12-21 05:51:45 START installation of SD Installer N 2015-12-21 05:51:45 Successful installation of SD Installer I 2015-12-21 05:51:45 END installation of SD Installer I 2015-12-21 05:51:45 START installation of Event Agent E 2015-12-21 05:53:28 Failed installation of Event Agent I 2015-12-21 05:53:28 END installation of Event Agent</pre> At the bottom, a message states: 'The System Management installation / registration / uninstallation may take several minutes. Please be patient while System Management is being installed / registered / uninstalled'. Navigation buttons 'Finish', 'Reset', and 'Go' are at the bottom right.	Functions	Installed	Status	SD- Asset- Management Agent Software Delivery - / Asset - Management	☒	not running	Event Management Agent Event Management	☒	not installed	Go Active	☐		Go Online	☐		☐
Functions	Installed	Status														
SD- Asset- Management Agent Software Delivery - / Asset - Management	☒	not running														
Event Management Agent Event Management	☒	not installed														
Go Active	☐															
Go Online	☐															
Click Go.	☐															

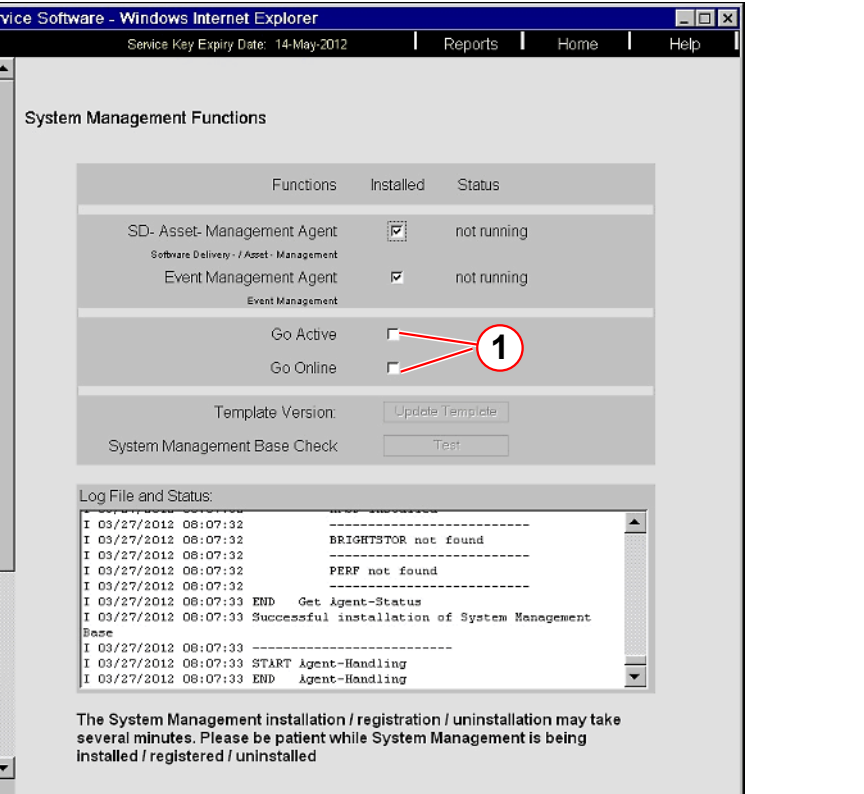
ACTION	MRAWP/MRWP
Now both Agents are installed and system management can be activated For the activation set the following configuration: Mark the checkboxes "Go Active" and "Go Online". (The other checkboxes have already been marked). <i>Fig. 178: System Management activation</i>  (1) Mark both checkboxes In the Log File and Status area, the current status of System Management activation/installation is displayed.	☐
Click Go.	☐
Note: The installation/registration will take approximately 10 to 15 minutes.	

ACTION

MRAWP/MRWP

- If activation is successful, the following is displayed:

Fig. 179: System Management activation OK

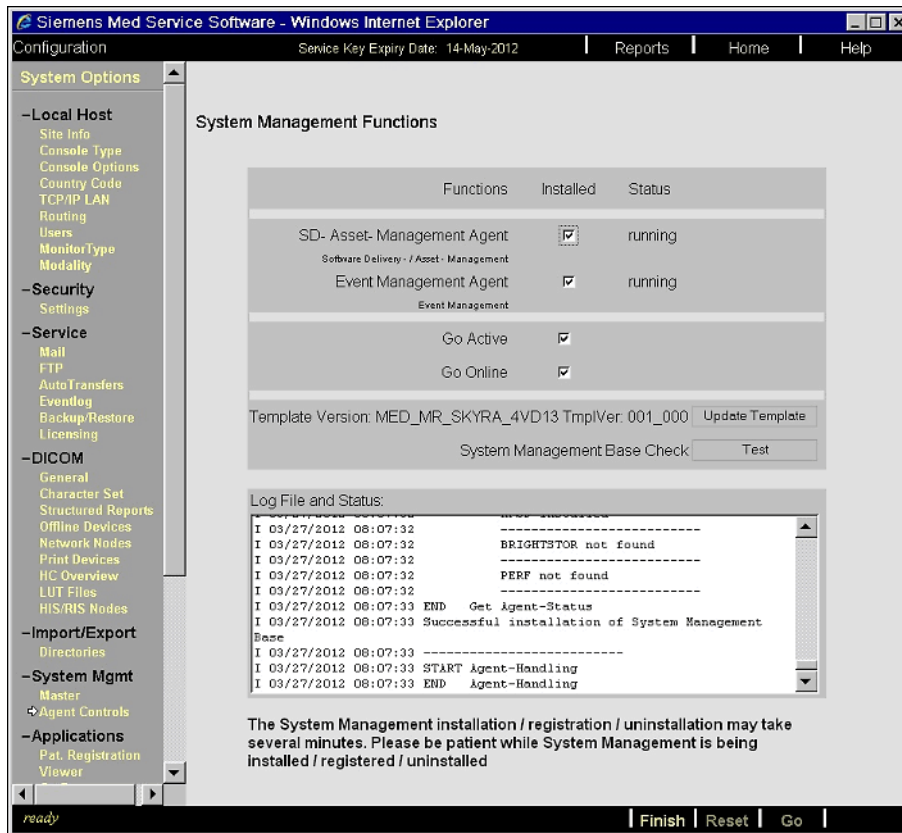


Select **Finish** to continue

- Click **Finish** and close the local service UI.


In the event of problems (installation/activation fails a second time), contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/Technical Support Group.

- Switch embedded control back into normal mode.

Fig. 180: MREmbeddedControlGUI update mode

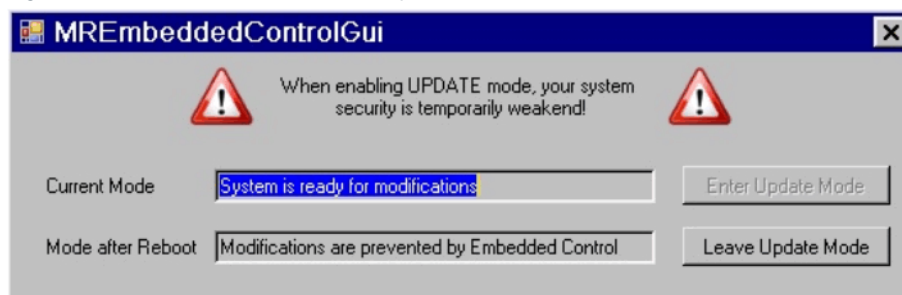




Press the button **Leave Update Mode** and confirm the displayed popup with **OK**

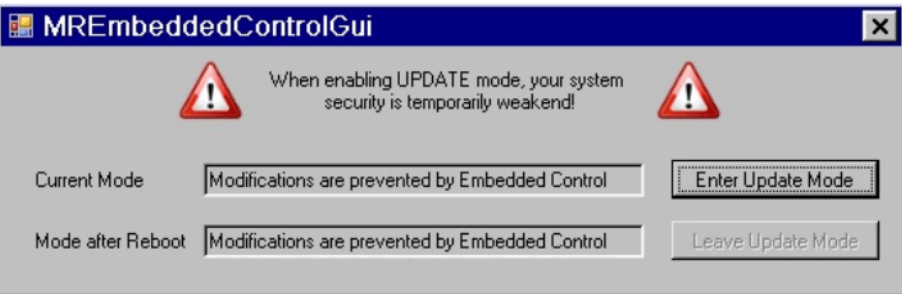ACTION	MRAWP/MRWP
Now EmbeddedControl is in normal mode. Fig. 181: MREmbeddedControlGUI 	☐
Close the MREmbeddedControlGUI by pressing the X	☐

4.1.10 Set Embedded Control in update mode



4.1.10.1 Logged in as user meduser

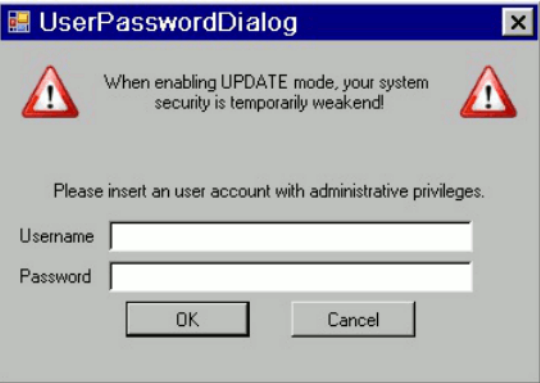ACTION	MRAWP/MRWP
Enable Advanced User	☐
Start the MREmbeddedControlGUI from the start menu Start ⇒ All Programs ⇒ MREmbeddedControl ⇒ MREmbeddedControlGUI	☐
In the displayed popup the user account with administrative privileges is requested. Fig. 182: EC User Password request  Type in as user the administrative account (e.g. syadmin) and the password.	☐

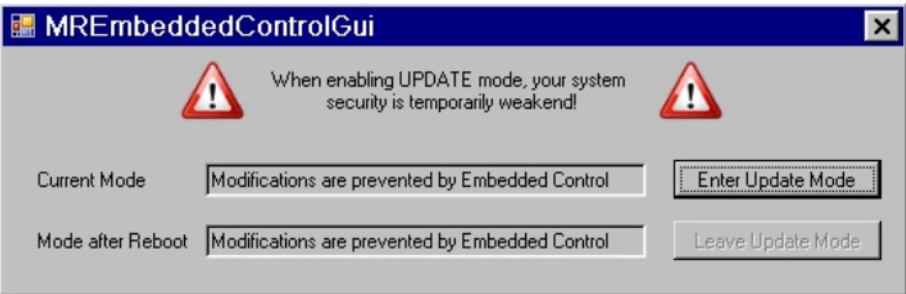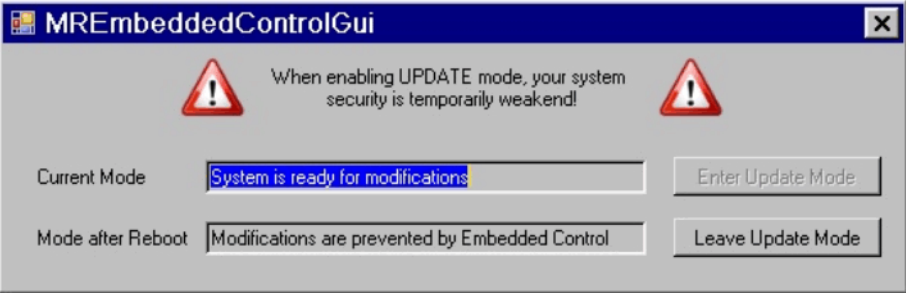ACTION	MRAWP/MRWP
Now the MREmbbedControlGUI is displayed <i>Fig. 183: MREmbbedControlGUI</i> 	☐
Press the button Enter Update Mode to switch embedded control into the update mode.	☐
Confirm the displayed popup with OK	☐
Now switch embedded control is in update mode. <i>Fig. 184: MREmbbedControlGUI update mode</i> 	☐
Switch embedded control back into normal mode, Press the button Leave Update Mode and confirm the displayed popup with OK	☐
Close the MREmbbedControlGUI by pressing the X	☐

4.1.10.2 Logged in as user syadmin (administrative privileges)

ACTION	MRAWP/MRWP
Enable Advanced User	☐
Start the MREmbbedControlGUI from the start menu Start ⇒ All Programs ⇒ MREmbbedControl ⇒ MREmbbedControlGUI	☐

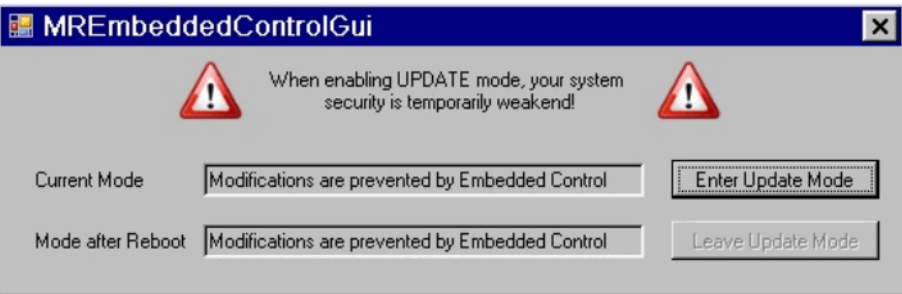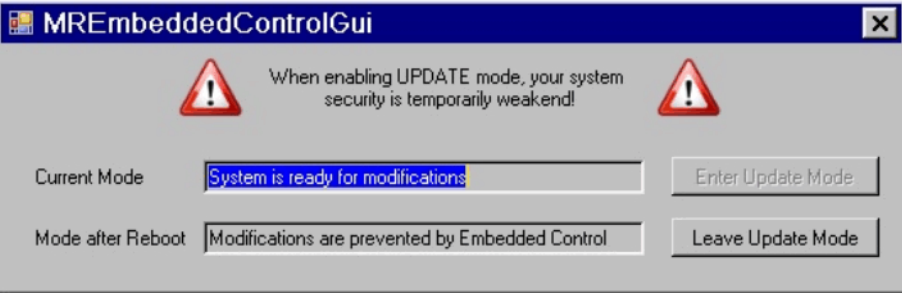ACTION	MRAWP/MRWP
Now the MREmbbedControlGUI is displayed Fig. 185: MREmbbedControlGUI 	☐
Press the button Enter Update Mode to switch embedded control into the update mode.	☐
Confirm the displayed popup with OK	☐
Now switch embedded control is in update mode. Fig. 186: MREmbbedControlGUI update mode 	☐
Switch embedded control back into normal mode, Press the button Leave Update Mode and confirm the displayed popup with OK	☐
Close the MREmbbedControlGUI by pressing the X	☐

4.1.11 Checking and setting the MOXA - Ports for HPz440



The tool "enable hidden dev" is required. Please download it in advance (→ enable hidden Dev, <https://productinformation.healthcare.siemens.com/FolderFile.ashx?filePath=XFRvb2xzXGVuYWJsZV9oaWRkZW5fRGV2LmV4ZQ%3d%3d&uuid=53b2f6bb-f1b0-4df9-bb77-dc5927732690>)





Embedded Control must be in update mode. Is this not the case it is not possible to extract the downloaded file and the execution of the script enable_hidden_Dev.bat is blocked. See for the update mode (→ Logged in as user syadmin (administrative privileges) / Page 147)




The Alarmbox serial cable has to be connected to the port of the MOXA USB device which is assigned to COM9 in the devicemanager.

The following procedure describes the correct settings of the com ports. This procedure has to be performed if the system standby function in the software (System Control) is not working.

ACTION	MRAWP
Ensure that you are logged in as with the administrative account (e.g. syadmin)" and embedde control is in update mode for it see (→ Logged in as user syadmin (administrative privileges) / Page 147)	☐
Insert the CD-ROM labeled "Tool enable hidden Dev" into the CD/DVD-ROM drive.	☐
- Start Windows Explorer: - Select START -> Programs -> Accessories -> Open Windows Explorer Windows Explorer starts.	☐
Navigate to the CD/DVD-ROM drive R:	☐
Extract the downloaded file enable_hidden_Dev.exe by double-clicking it.	☐
One file is extracted to C:\Temp\tools directory	☐
Invoke a command shell from the start menu (Start All ⇒ Programs Accessories ⇒ Command Prompt)	☐
Execute the enable_hidden_Dev.bat script by typing in the cmd tool window: C:\Temp\tools\enable_hidden_Dev.bat ↵	☐

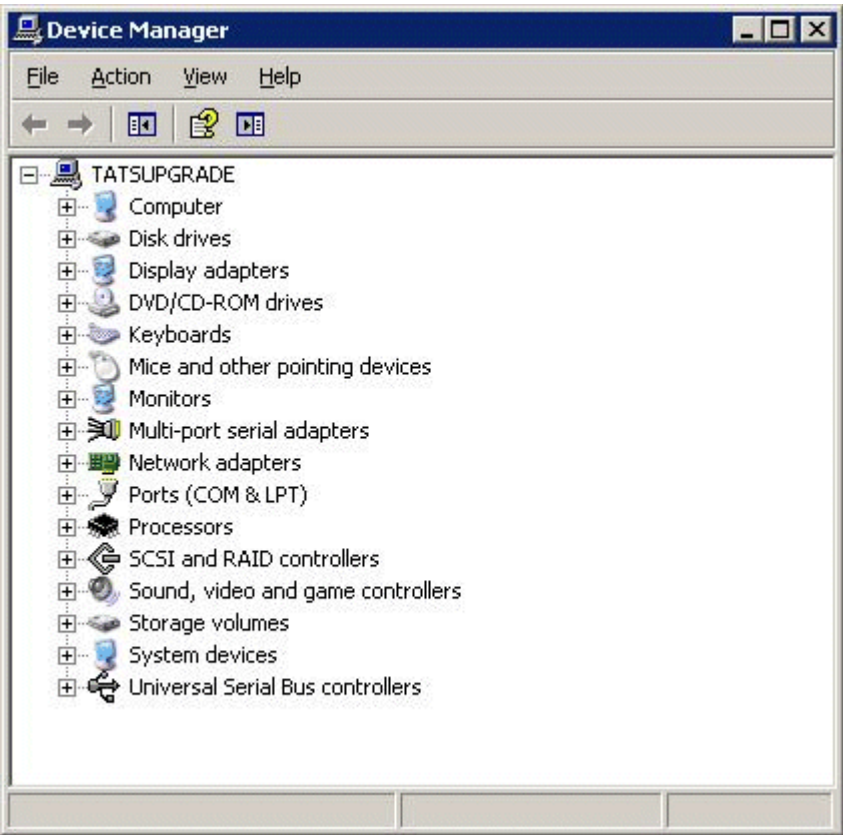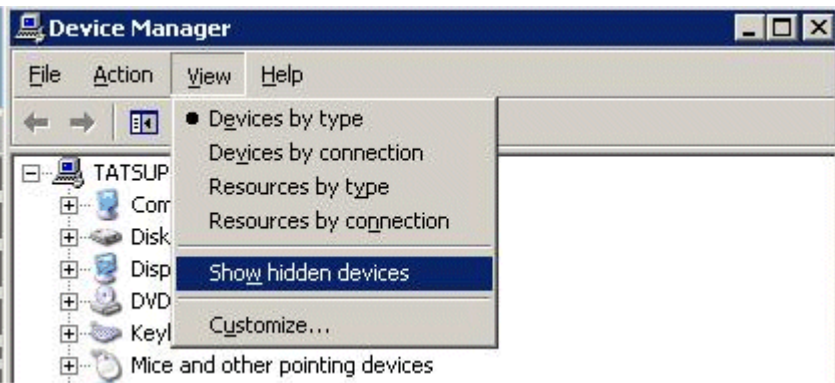ACTION	MRAWP
This displays the Device Manager. Fig. 187: Device manager A screenshot of the Windows Device Manager window. The title bar says "Device Manager". The menu bar has "File", "Action", "View", and "Help". Below the menu bar is a toolbar with icons for back, forward, and a help icon. The main area shows a tree view of hardware categories under the name "TATSUPGRADE". The categories are: Computer, Disk drives, Display adapters, DVD/CD-ROM drives, Keyboards, Mice and other pointing devices, Monitors, Multi-port serial adapters, Network adapters, Ports (COM & LPT), Processors, SCSI and RAID controllers, Sound, video and game controllers, Storage volumes, System devices, and Universal Serial Bus controllers. Each category has a plus sign icon next to it, indicating it can be expanded.	☐
Enable the view that shows hidden devices. Select View -> Show hidden devices Fig. 188: Device manager - show hidden devices A screenshot of the Windows Device Manager window with the "View" menu open. The menu options are: "Devices by type", "Devices by connection", "Resources by type", "Resources by connection", "Show hidden devices" (which is highlighted with a blue background), and "Customize...". The tree view on the left shows the same categories as in Fig. 187, but only the first few are visible: "TATSUP", "Com", "Disk", "Disp", "DVD", "Key", and "Mice and other pointing devices".	☐
- Expand the entry Multi-port serial adapters in the device manager by clicking +. If you should get two entries of: - MOXA UPort 2210 Series Both entries for MOXA UPort 2210 Series must be removed.	☐

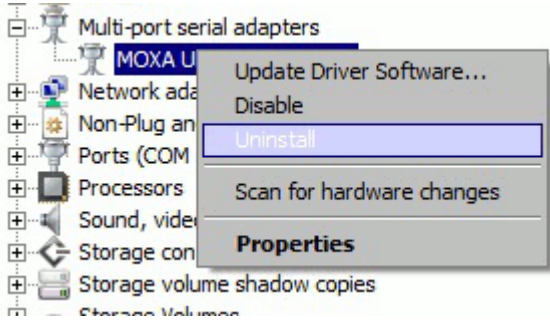
ACTION	MRAWP
■ To do this, proceed as follows: - Select the first entry MOXA UPort 2210 Series and press the right mouse button. <i>Fig. 189: Uninstall device</i>  - Select Uninstall in the displayed menu. - Confirm the uninstallation request pop-up with OK. This will uninstall the device. - Select the next entry for MOXA Uport 2210 Series and press the right mouse button. - Select the Uninstall entry in the displayed menu. - Confirm the uninstallation request pop-up with OK. This will uninstall the device.	☐
■ Reinstall the network devices. To do so, select the entry Action -> Scan for hardware changes in the device manager menu. <i>Fig. 190: Device manager - reinstall devices</i> 	☐

ACTION	MRAWP
■ The network devices are now reinstalled. <i>Fig. 191: Device manager - Multi-port serial adapter reinstalled</i>  The screenshot shows the Windows Device Manager window. The title bar reads 'Device Manager'. The menu bar includes 'File', 'Action', 'View', and 'Help'. Below the menu is a toolbar with icons for back, forward, refresh, help, and search. The main pane shows a tree view of device categories. The 'Multi-port serial adapters' category is expanded and highlighted with a blue selection box. Under this category, the 'MOXA UPort 2210 Series' is listed. Other categories visible include Batteries, Computer, Disk drives, Display adapters, DVD/CD-ROM drives, Human Interface Devices, Keyboards, Mice and other pointing devices, Monitors, Network adapters, Ports (COM & LPT), Processors, Sound, video and game controllers, Storage controllers, System devices, and Universal Serial Bus controllers.	☐

ACTION

- Right-click "MOXA UPort...." and select "Properties"

Fig. 192: Device manager - MOXA port check

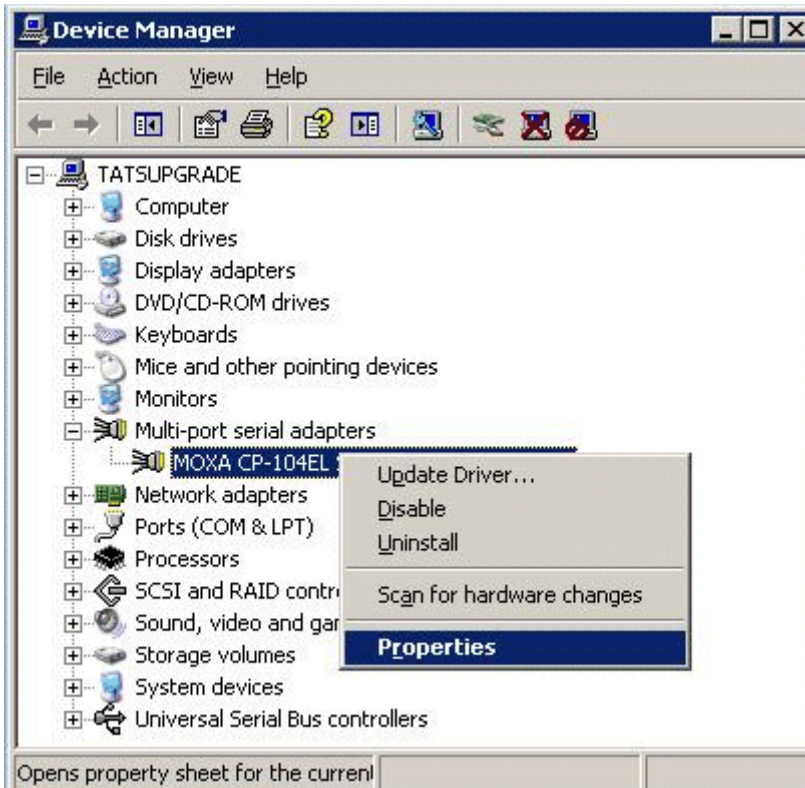




MRAWP



MAGNETOM	M6-000.896.65.01.02 12.18	Page 154 of 162	© Siemens Healthineers, 2018 Restricted
----------	------------------------------	-----------------	--

ACTION	MRAWP
■ Select Port 1 and Click "Port setting" - The "Auto Enumerating COM Number" checkbox is marked and so are the other checkboxes. - Set the Port Number of Port 1 to COM 6 Fig. 194: MOXA port 1 configuration Port 1 Port NumberCOM6 (current) ☒ Auto enumerating COM number Friendly NameMOXA UPort COM ☒ Set the change to all ports Fast FlushEnable ☒ Set the change to all ports InterfaceRS-232 ☒ Set the change to all ports OK Cancel Press "OK"	☐

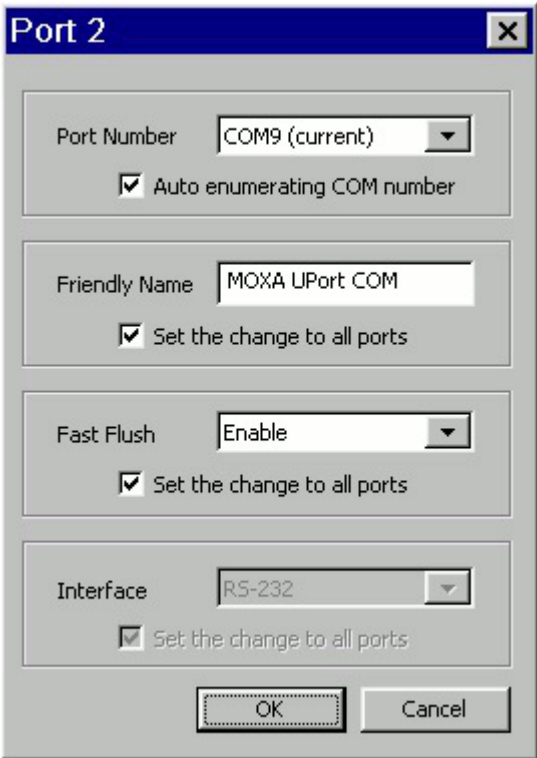ACTION	MRAWP
■ Select Port 2 and Click "Port setting" - The "Auto Enumerating COM Number" checkbox is marked and so are the other checkboxes. - Set the Port Number of Port 2 to COM 9 Fig. 195: MOXA port 2 configuration  Press "OK"	☐

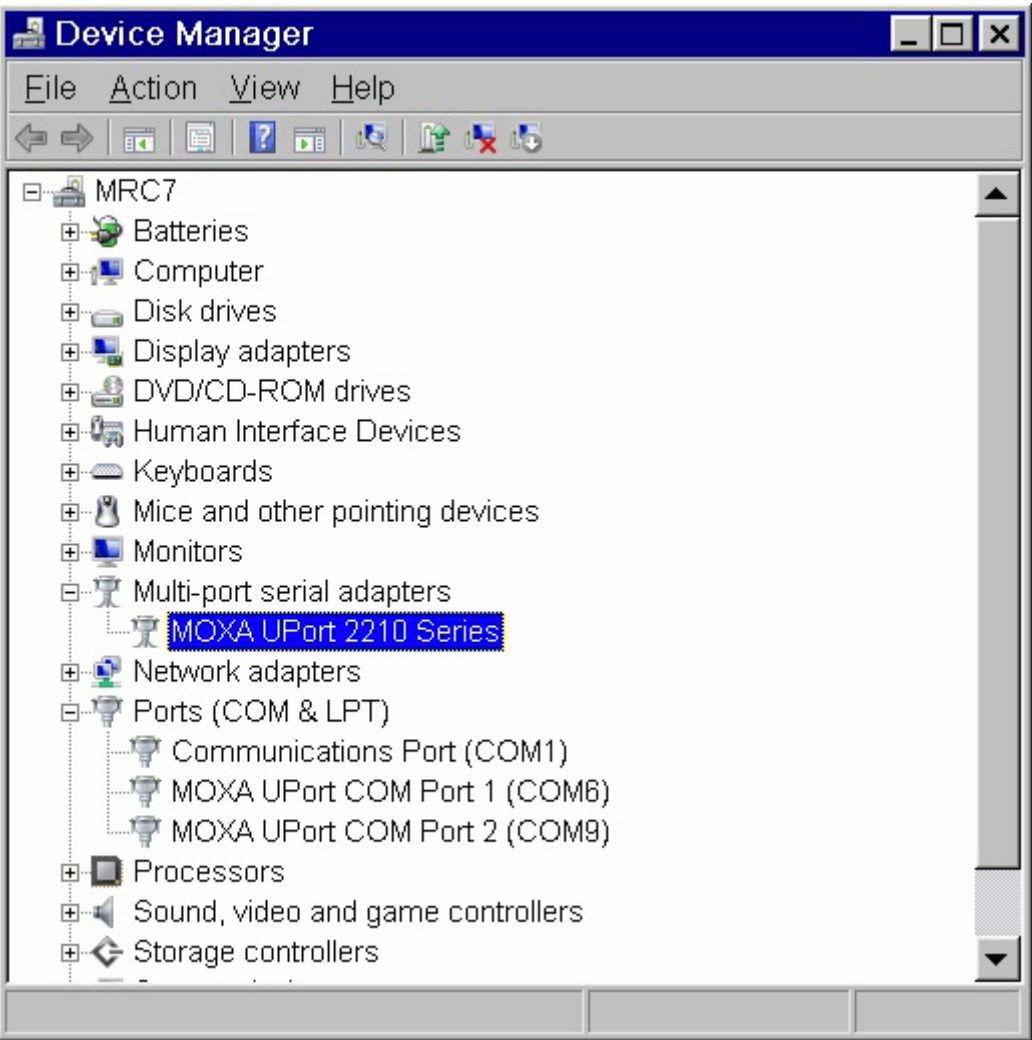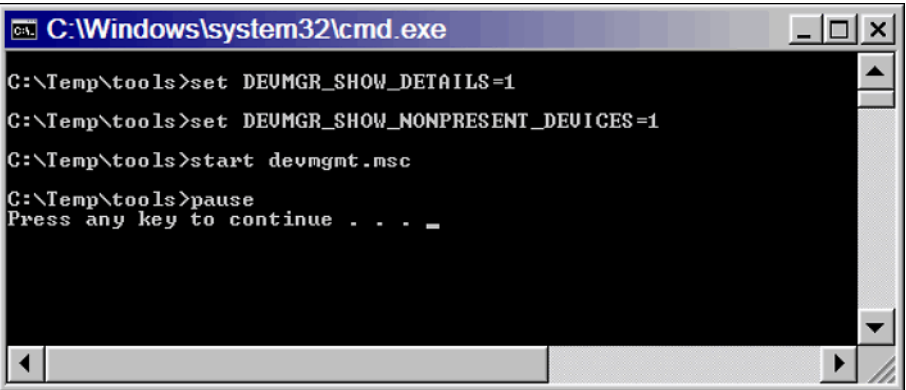ACTION	MRAWP
You can now check to ensure that these settings under "Ports (COM & LPT)" are automatically set to COM 6 and COM 9 as shown in the picture below. Fig. 196: Device manager - MOXA COM port assignment 	☐
Close the device manager.	☐
Close the following cmd tool by pressing the space bar. Fig. 197: CMD window 	☐

ACTION	MRAWP
■ Reboot the system.	☐
■ Login as user meduser if autologin is disabled and wait until the syngo main UI is up.	☐
■ Perform the Standby function check again.	☐

5.1 VE11D Evolve

5.1.1 Reason for the Update

This Update Instruction is for obligating Evolve update. It provides the compatible hardware that is required (host(s) and MaRS) and the software syngo MR E11D according to the customer's Evolve agreement. The software VE11D does contain:

- Significant enhancement in security configuration for secure scanner operation in hospital network.
- syngo Security Enhanced

In case the previous SW version was VD13 the VE11D provides additionally the new Dot Cockpit for easy Dot Engine configuration together with several workflow and performance enhancements:

- DotGo
- Dot Cockpit

With the available options:

- StarVibe (FREEZEit)
- TWIST-Vibe (FREEZEit)
- LiverLab
- MyoMaps
- Dot Engines (Spine, Large Joint, Breast)

5.1.2 Protocol

This update has been completed successfully. The customer has been informed of the benefits or improvements resulting from this update.

Customer:

Customer No.: Func. Location

Material No. Serial number
of the system: of the system:

Notification:

Remark:

Remark:

Country: Site:

Performed by: Telephone:

Date: Signature:



File the protocol in the corresponding binder.

There are no Hazard IDs in this document.

- Restricted - All documents may only be used by authorized personnel for rendering services on Siemens Healthcare Products. Any document in electronic form may be printed once. Copy and distribution of electronic documents and hardcopies is prohibited. Offenders will be liable for damages. All other rights are reserved.

siemens-healthineers.com

Siemens Healthineers Headquarters
Siemens Healthcare GmbH
Henkestr. 127
91052 Erlangen
Germany
Telephone: +49 9131 84-0
siemens-healthineers.com

Print No.: M6-000.896.65.01.02 | Replaces: n.a.
Doc. Gen. Date: 12.18 | Language: English
© Siemens Healthineers, 2018