



Power Supply Replacement Instructions

GEHC FRU 5390239

1. PROCEDURE

This procedure is applicable for replacing the 12Vdc power supply in the following panel:
GEHC P/N 5161287 Discovery MR750 Main Disconnect Panel which may also be referred to as the DVMR Main Disconnect Panel.

Reminder:



RISK OF ELECTRICAL HAZARD

Ensure all power sources feeding this PDB panel have been turned off and securely locked-out and tagged-out using the GEHC and hospital Lock-Out/Tag-Out (LOTO) procedures.

Depending on local rules, use protective equipment and tools. (Gloves and Arc Flash mask)

Required Tools & Equipment

- Screwdrivers, various sizes, slotted & Phillips
- Digital Voltmeter

2. POWER SUPPLY INSTALLATION



- Ensure the main disconnect is turned off and properly Locked-Out and Tagged-Out.
- Open the enclosure door.
- Open CB5, removing the load circuit from the power supply.
- Locate the power supply.

The actual power supply may be mounted in one of two configurations:

It may be mounted to DIN rail, fastened directly to the enclosure subpanel.

OR

It may be mounted to DIN rail on a small adapter plate attached to the top of the AC duplex receptacle.

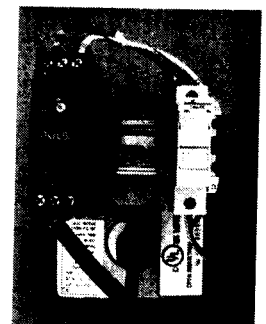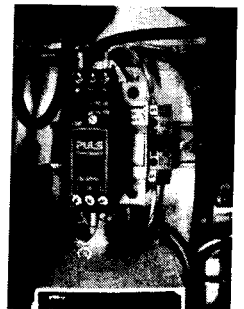



- e. Disconnect the wires from the power supply line side.

DIN Rail Mounted P. S.	
P.S Terminal	Wire # - Color
L	3 - Red
N	2 - White
GND	- Grn/Yel

Outlet Mounted P.S.	
P.S Terminal	Wire # - Color
L	Power Cord - Brown
N	Power Cord - Blue
GND	Power Cord - Grn/Yel

- f. Disconnect the wires from the power supply DC output side.

P.S Terminal	Wire # - Color
+	11A - Blue
-	12 - White/Blue

- g. Remove the old power supply from the DIN rail. Release the power supply retaining clip with a small screwdriver and rotate the power supply up and away from the DIN rail.
- h. Install the new power supply on the DIN rail.
- i. Engage power supply on to the top side of the DIN rail.
- j. Rotate the power supply downward, snapping the power supply onto the bottom side of the DIN rail.
- k. Re-connect the wires to the new power supply according to the wire charts above.



3. VERIFY THE POWER SUPPLY OUTPUT VOLTAGE

- a. Ensure the Load equipment, the PDU and Heat Exchange Cabinet are turned off.
- b. Remove the Lock-out/Tag-out devices.
- c. Apply the 480 Vac source power to the panel. Turn the Main Circuit Breaker, CB1, ON. The green indicator light on the power supply should be illuminated.
- d. Ensure the FU4 fuseholder is still open.
- e. While measuring the DC voltage between the load side of the power supply [wires #11A (+) wire #12, (-), adjust the power supply output voltage for 13.45 to 13.55 VDC.
- f. Close the FU4 fuseholder.
- g. Close the enclosure door and place the system back into operation.
- h. Dispose of the old power supply using proper GEHC disposal procedures.

Rev	Description of Change	Changed By	Approved By	Date
A	Initial Release	JFH		