

MAGNETOM

Replacement of Parts System

Patient Handling

Avanto
Espree

Document Version

Siemens reserves the right to change its products and services at any time.

In addition, manuals are subject to change without notice. The hardcopy documents correspond to the version at the time of system delivery and/or printout. Versions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order a current version or refer to our website <http://www.healthcare.siemens.com>.

Disclaimer

Siemens provides this documentation “as is” without the assumption of any liability under any theory of law.

The content described herein requires superior understanding of our equipment and may only be performed by qualified personnel who are specially trained for such installation and/or service.

Copyright

“© Siemens, 2006” refers to the copyright of a Siemens entity such as:

Siemens Healthcare GmbH - Germany

Siemens Aktiengesellschaft - Germany

Siemens Shenzhen Magnetic Resonance Ltd. - China

Siemens Shanghai Medical Equipment Ltd. - China

Siemens Healthcare Private Ltd. - India

Siemens Medical Solutions USA Inc. - USA

Siemens Healthcare Diagnostics Inc. - USA and/or

Siemens Healthcare Diagnostics Products GmbH - Germany

1	General	8
2	FRU Overview	10
3	Patient table	16
3.1	Guide rails	16
3.1.1	General	16
3.1.2	Tools and Time required	16
3.1.3	Preliminary	17
3.1.4	Disassembly	17
3.1.5	Assembly	27
3.1.6	Concluding steps	33
3.1.7	Adjustments	33
3.1.8	Final checks	33
3.2	Lifting motor / toothed belt of vertical drive	34
3.2.1	Tools and time required	34
3.2.2	Preliminary	34
3.2.3	Disassembly	34
3.2.4	Assembly	36
3.2.5	Concluding steps	36
3.2.6	Adjustments	36
3.2.7	Final checks	36
3.3	Spindle	37
3.3.1	General	37
3.3.2	Spare part items of spindle	39
3.3.3	Principle of spindle (old version)	41
3.3.4	Tools and time required	41
3.3.5	Preliminary	42
3.3.6	Disassembly method I with table transport device	42
3.3.7	Disassembly method II, with a massive wood beam	51
3.3.8	Assembly method I	61
3.3.9	Assembly method II	64
3.3.10	Concluding steps	66
3.3.11	Adjustments	66
3.3.12	Final checks	66
3.4	Brake disk	67
3.4.1	General	67
3.4.2	Tools and time required	69
3.4.3	Preliminary	69
3.4.4	Disassembly	69
3.4.5	Assembly	78
3.4.6	Concluding steps	78
3.4.7	Final checks	78
3.5	Tabletop	79
3.5.1	General	79
3.5.2	Tools and time required	79
3.5.3	Preliminary	79
3.5.4	Disassembly	80
3.5.5	Assembly	84

3.5.6	Concluding steps	85
3.5.7	Adjustments	86
3.5.8	Final checks	86
3.6	Body of the energy chain	87
3.6.1	General	87
3.6.2	Tools and time required	87
3.6.3	Preliminary	87
3.6.4	Disassembly	87
3.6.5	Assembly	91
3.6.6	Concluding steps	92
3.6.7	Adjustments	92
3.6.8	Final checks	92
3.7	Guide rollers	93
3.7.1	General	93
3.7.2	Tools and time required	93
3.7.3	Preliminary	93
3.7.4	Disassembly	93
3.7.5	Assembly	96
3.7.6	Concluding steps	96
3.7.7	Adjustments	96
3.7.8	Final checks	96
3.8	MDSD device	97
3.8.1	General	97
3.8.2	Tools and time required	97
3.8.3	Preliminary	98
3.8.4	Disassembly	98
3.8.5	Assembly	116
3.8.6	Adjustments	124
3.8.7	Concluding steps	124
3.8.8	Adjustments	125
3.8.9	Final checks	125
3.9	Horizontal drive unit	126
3.9.1	General	126
3.9.2	Tools and time required	126
3.9.3	Preliminary	126
3.9.4	Disassembly	126
3.9.5	Assembly	146
3.9.6	Concluding steps	155
3.9.7	Adjustments, not valid for Trio A Tim Systems	156
3.9.8	Final checks	157
3.10	Toothed belt of horizontal drive	158
3.10.1	General	158
3.10.2	Tools and time required	158
3.10.3	Preliminary	158
3.10.4	Disassembly	158
3.10.5	Assembly	167
3.10.6	Adjustments	170
3.10.7	Concluding steps	171
3.10.8	Final checks	171
3.11	Sound transducer	172

3.11.1	General	172
3.11.2	Tools and time required	172
3.11.3	Preliminary	172
3.11.4	Disassembly	173
3.11.5	Assembly	178
3.11.6	Concluding steps	179
3.11.7	Adjustments	179
3.11.8	Final checks	179
3.12	RF coil plugs	180
3.12.1	General	180
3.12.2	Tools and time required	180
3.12.3	Preliminary	180
3.12.4	Disassembly	180
3.12.5	Assembly	184
3.12.6	Concluding steps	184
3.12.7	Final checks	185
3.12.8	Disassembly	185
3.12.9	Assembly	191
3.12.10	Concluding steps	193
3.12.11	Final checks	193
3.13	Control switch S30	194
3.13.1	General	194
3.13.2	Tools and time required	194
3.13.3	Preliminary	194
3.13.4	Disassembly	196
3.13.5	Assembly	201
3.13.6	Concluding steps	202
3.13.7	Final checks	202
3.14	Cable replacement	204
3.14.1	General	204
3.14.2	Tools and time required	204
3.14.3	Preliminary	204
3.14.4	Disassembly	204
3.15	Control panel	216
3.15.1	General	216
3.15.2	Tools and time required	216
3.15.3	Preliminary	216
3.15.4	Disassembly	216
3.15.5	Assembly	225
3.15.6	Concluding steps	225
3.15.7	Adjustments	225
3.15.8	Final checks	225
3.16	Patient light Espree	226
3.16.1	General	226
3.16.2	Tools and time required	226
3.16.3	Preliminary	226
3.16.4	Disassembly	227
3.16.5	Assembly	229
3.16.6	Concluding steps	229
3.16.7	Adjustments	229

3.16.8	Final checks	229
3.17	Patient light Avanto	230
3.17.1	General	230
3.17.2	Tools and time required	230
3.17.3	Disassembly	230
3.17.4	Assembly	231
3.17.5	Concluding steps	231
3.17.6	Adjustments	231
3.17.7	Final checks	231
3.18	Light localizer / PMU antenna / microphone	232
3.18.1	General	232
3.18.2	Tools and time required	232
3.18.3	Preliminary	232
3.18.4	Disassembly	233
3.18.5	Assembly	236
3.18.6	Concluding steps	236
3.18.7	Adjustments	236
3.18.8	Final checks	237
4	Trolley	238
4.1	Locking Hooks (Mechanical)	238
4.1.1	General	238
4.1.2	Tools and Time required	238
4.1.3	Preliminary	238
4.1.4	Disassembly	239
4.1.5	Assembly	246
4.1.6	Final Steps	248
4.1.7	Final Checks	248
4.2	Keypads / PCBs	249
4.2.1	General	249
4.2.2	Tools and time required	249
4.2.3	Preliminary	249
4.2.4	Disassembly	249
4.2.5	Assembly	251
4.2.6	Concluding steps	251
4.2.7	Final checks	252
4.3	Wheel Reinforcement Kit	253
4.3.1	General	253
4.3.2	Tools and Time required	253
4.3.3	Preliminary	253
4.3.4	Installing the Wheel Reinforcement Kit	254
4.3.5	Final Checks	262
5	PMU	263
5.1	PDAU	263
5.1.1	General	263
5.1.2	Tools and time required	263
5.1.3	Preliminary steps	263
5.1.4	Disassembly	263
5.1.5	Assembly	265

5.1.6	Concluding steps	265
5.1.7	Adjustments	265
5.1.8	Final checks	265
6	Video monitoring / intercom	266
6.1	Loudspeaker/Microphone	266
6.1.1	General	266
6.1.2	Tools and time required	266
6.1.3	Preliminary	266
6.1.4	Disassembly	266
6.1.5	Assembly	267
6.1.6	Concluding steps	267
6.1.7	Adjustments	268
6.1.8	Final checks	268
6.1.9	Disassembly	268
6.1.10	Assembly	270
6.1.11	Concluding steps	270
6.1.12	Adjustments	270
6.1.13	Final checks	270
6.2	CCD camera	271
6.2.1	General	271
6.2.2	Tools and time required	271
6.2.3	Preliminary	271
6.2.4	Disassembly	271
6.2.5	Concluding steps	273
6.2.6	Adjustments	273
6.2.7	Final checks	273
7	Changes to previous version	274
8	List of Hazard IDs	275

This chapter describes the replacement procedure for all patient handling components of the MR system and consists of the following:

- Patient table, **ptab**
- Removable tabletop (option), **trolley**
- Physiological measurement unit **PMU**
- Video monitoring and intercom



CAUTION

Strong magnetic field!

Magnetic objects may cause injury.

- » Ramp down the magnet or maintain a safe distance to the magnetic field.



CAUTION

Strong magnetic field!

Magnetic objects may cause injury.

- » Hold the ferromagnetic part with both hands while moving it to the correct position and to remove the unit directly out of the field!



Switching **off** the patient table via circuit breaker F22 may cause the patient table to re-load.



For disassembling/assembling the patient table covers, we recommend using a midsize Philips screwdriver



Please consider the different length of screws during assembly of the patient table covers



If any end stop of the patient table is removed due to disassembly, the end stop must be reattached to the original position.



Adhere to ESD regulations and use the appropriate tool set: Anti-static wristband and grounded pad. (TDK equipment Type 8005 from 3M Deutschland GmbH)



Protective ground connection test.

If a protective earth connection has to be disconnected for replacing a component, measure the protective earth connection resistance $\leq 100 \text{ m}\Omega$ with a ground test meter (Fluke 1503, or equivalent) between the protective earth connector on the "sub-component" (e.g. direct position drive) and the protective earth wire connected to the "main component" (e.g. patient table).

Document the successful measurement in the service report (measured value and measured equipment) as required by IEC 62353.



CAUTION

Possible damage of the spindle nut.

- » When using the table transport device or a wooden support beam, use the manual crank only until the load is removed from the spindle nut. Do not continue to crank down manually without a load to the spindle nut.

This chapter described all FRUs in the Patient Handling System that can be replaced.

Patient table					
FRU	Time	Instructions	Adjustments	Final checks	Tool
Guide rails	120 min.	(→ Patient table / Page 16)	■ n.a.	■ Vertical table movement	Non-magnetic tool kit Electric drill Hammer Plastic block Table transport device Loctite 221
Lifting motor/vertical	30 min.	(→ Lifting motor / toothed belt of vertical drive / Page 34)	■ n.a.	■ Vertical table movement	Non-magnetic tool kit
Toothed belt - vertical drive	30 min.	(→ Lifting motor / toothed belt of vertical drive / Page 34)	■ n.a.	■ Vertical table movement	Non-magnetic tool kit
Spindle	120 min./60 min.	(→ Spindle / Page 37)	■ n.a.	■ Vertical table movement	Non-magnetic tool kit Table transport device Loctite 221
Brake Disk	60 min.	(→ Brake disk / Page 67)	■ n.a.	■ Vertical table movement	Non-magnetic tool kit Loctite 221
Tabletop plate	60 min.	(→ Tabletop / Page 79)	■ n.a.	■ Horizontal table movement	Non-magnetic tool kit
Energy chain	30/45 min.	(→ Body of the energy chain / Page 87)	■ n.a.	■ Horizontal table movement	Non-magnetic tool kit
Guide rollers	30/45 min.	(→ Guide rollers / Page 93)	■ n.a.	■ Horizontal table movement	Non-magnetic tool kit Loctite 221

Patient table					
FRU	Time	Instructions	Adjustments	Final checks	Tool
MDSD device	120 min.	(→ MDSD device / Page 97)	■ n.a.	■ Horizontal table movement	Non-magnetic tool kit Non-magnetic torque wrench Flat open-end spanner, 22/27 mm Flat blade screwdriver 4 mm Philips screwdriver, medium size
Horizontal drive unit	120 min.	(→ Horizontal drive unit / Page 126)	■ n.a.	■ Horizontal table movement	Non-magnetic tool kit Non-magnetic torque wrench Flat open-end spanner, 22 mm Flat blade screwdriver 4 mm Philips screwdriver, medium size Non-magnetic torque wrench Loctite 221
Toothed belt of horizontal drive	60 min.	(→ Toothed belt of horizontal drive / Page 158)	■ n.a.	■ Horizontal table movement	Non-magnetic tool kit Philips screwdriver, medium size
Sound transducer	60 min.	(→ Sound transducer / Page 172)	■ n.a.	■ Table functions, patient communication	Non-magnetic tool kit Philips screwdriver, medium size
Circuit boards of patient table electronics, A4120, A4140, A4110	30 to 90 min.	Check status of DIP switches	■ n.a.	■ Table functions	Non-magnetic tool kit

Patient table					
FRU	Time	Instructions	Adjustments	Final checks	Tool
RF coil plugs/cap kit	30 min. for each side	(→ RF coil plugs / Page 180)	■ n.a.	■ Quality Assurance / Coil check	Non-magnetic tool kit Service plug
Switch S30	60 min.	(→ Control switch S30 / Page 194)	■ n.a.	■ Switch test in Service mode	Non-magnetic tool kit
Cable	180 min.	(→ Cable replacement / Page 204)	■ n.a.	■ Function test of replaced cable	Non-magnetic tool kit
Transformer A4114	30 min.+ Magnet ramp time	■ n.a.	■ n.a.	■ Function test of the transformer - Check table movement/ horizontal and vertical	Caution: Magnetic part. Ramp down the magnet.
Control panel Espre	20 min.	(→ Control panel / Page 216)	■ n.a.	■ Function test of control panel	Non-magnetic tool kit
Keypad foil Espre	25 min.	(→ Control panel / Page 216)	■ n.a.	■ Function test of control panel	Non-magnetic tool kit
Patient bore light Espre	30 min.	(→ Patient light Espre / Page 226)	■ n.a.	■ Function test of patient bore light	Non-magnetic tool kit
Patient light Avanto	10 min.	(→ Patient light Avanto / Page 230)	■ n.a.	■ Function test of patient bore light	Non-magnetic tool kit
Light localizer Espre	30 min.	(→ Light localizer / PMU antenna / microphone / Page 232)	■ Light localizer	■ Function test of light localizer	Non-magnetic tool kit

Trolley					
FRU	Time	Instructions	Adjustments	Final checks	Tool
Locking hooks	120 min.	(→ Locking Hooks (Mechanical) / Page 238)	■ n.a.	■ Check for smooth and safe trolley operation	Non-magnetic tool kit
Keypads/PCB	20 min.	(→ Keypads / PCBs / Page 249)	■ n.a.	■ Check smooth and save trolley operation	Non-magnetic tool kit

PMU					
FRU	Time	Instructions	Adjustments	Final checks	Tool
PDAU	20 min.	(→ PMU / Page 263)	■ n.a.	■ PMU Test	Non-magnetic tool kit
PERU	2 min.	n.a.	■ n.a.	■ PMU Test	n.a.
PPU	2 min.	n.a.	■ n.a.	■ PMU Test	n.a.
PMU antenna	25 min.	(→ Light localizer / PMU antenna / microphone / Page 232)	■ n.a.	■ PMU Test	n.a.
Sensor charger	5 min.	n.a.	■ n.a.	■ PMU Test	n.a.
Cable kit PMU Display	15 min.	n.a.	■ n.a.	■ PMU Test	n.a.

Video monitoring / intercom					
FRU	Time	Instructions	Adjustments	Final checks	Tool
Loudspeaker, Microphone Avanto	20 min.	(→ Video monitoring / intercom / Page 266)	■ n.a.	■ Acoustic test	Non-magnetic tool kit

Video monitoring / intercom					
FRU	Time	Instructions	Adjustments	Final checks	Tool
Microphone Espre	30 min.	(→ Video monitoring / intercom / Page 266)	■ n.a.	■ Acoustic test	Non-magnetic tool kit
Fan Espre - contains magnetic objects, maintain safe distance	30 min.	n.a.	■ n.a.	■ Functional test	Non-magnetic tool kit
Intercom central unit	10 min.	n.a.	■ n.a.	■ Bi-directional acoustic test, table stop button	Non-magnetic tool kit
Control panel intercom	5 min.	n.a.	■ n.a.	■ Bi-directional acoustic test, table stop button	Non-magnetic tool kit
Patient monitor	5 min.	n.a.	■ n.a.	■ Visual test	Non-magnetic tool kit
CCD Camera Avanto	15 min.	(→ CCD camera / Page 271)	■ n.a.	■ Visual test	Non-magnetic tool kit
Video Camera Espre	15 min.	(→ Start-up / MR06-010.814.02)	■ Optical adjustment/Start-up	■ Visual test	Non-magnetic tool kit

3.1 Guide rails

3.1.1 General

The table is moved in a vertical direction via four preassembled bearing systems above two guide rails. The guide rails and the bearings must be replaced together.

3.1.2 Tools and Time required

3.1.2.1 Time

The estimated time does **not include** ramping the magnet.
Approx. 120 min

3.1.2.2 Tools

- Loctite 221
- Non-magnetic tool kit

Fig. 1: Allen key



- 5 mm Allen key delivered with the guide rail kit.

- Table transport device.

Fig. 2: Table transport device


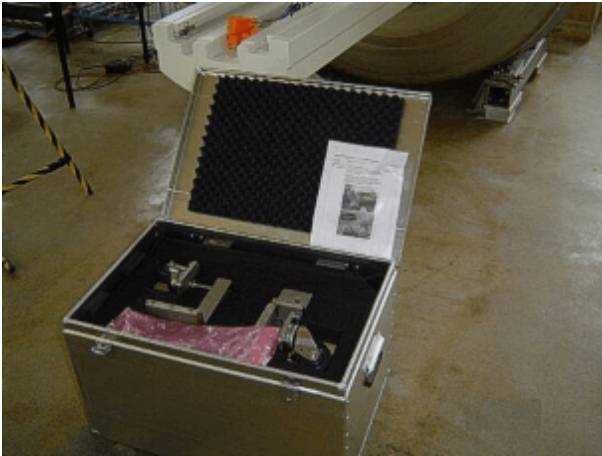


Fig. 3: Table transport device

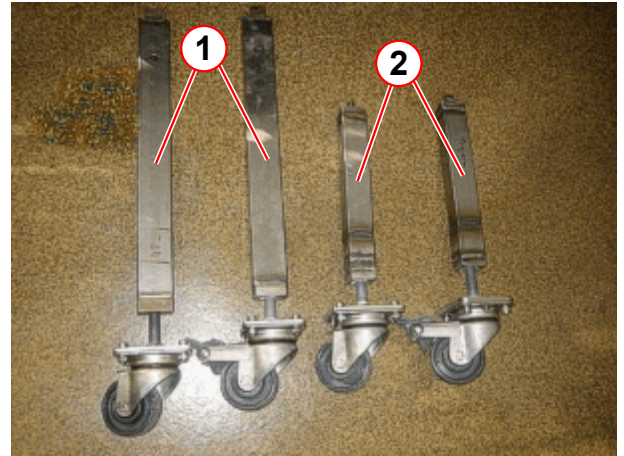

- (1) Transport roller 1
(2) Transport roller 2

3.1.3 Preliminary

- Ramp down the magnet.
- Switch OFF F22 (RF-room) in the ACC.



CAUTION

Strong magnetic field!

Magnetic objects may cause injury.

» Ramp down the magnet

3.1.4 Disassembly

3.1.4.1 Covers



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

- Remove the service access cover and the lifting column cover.


The removal of the service access cover and the lifting column cover is described in the REP/system, refer to (→ System / M6-020.841.10).

Fig. 4: Lower magnet, front cover



- Remove lower front cover of the magnet.


Fig. 5: Table, front cover

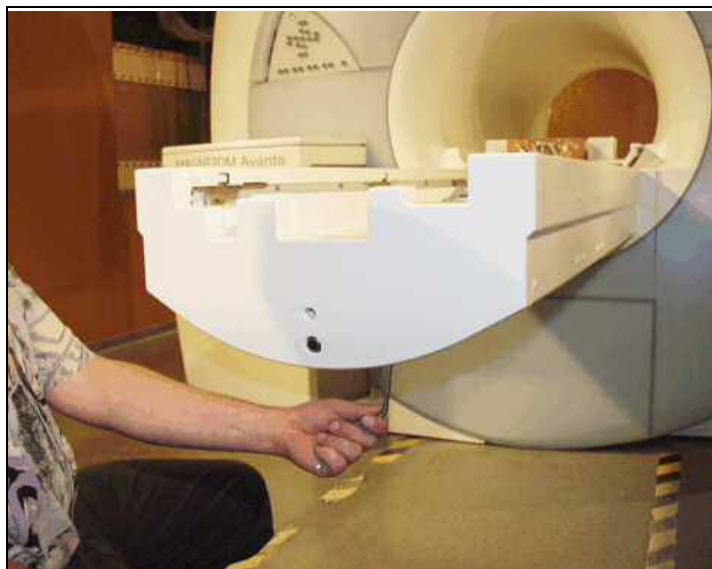

- Pull the table out of the magnet and remove the front cover of the table.

3.1.4.2 Bellows guide

Fig. 6: Table down, Allen screws

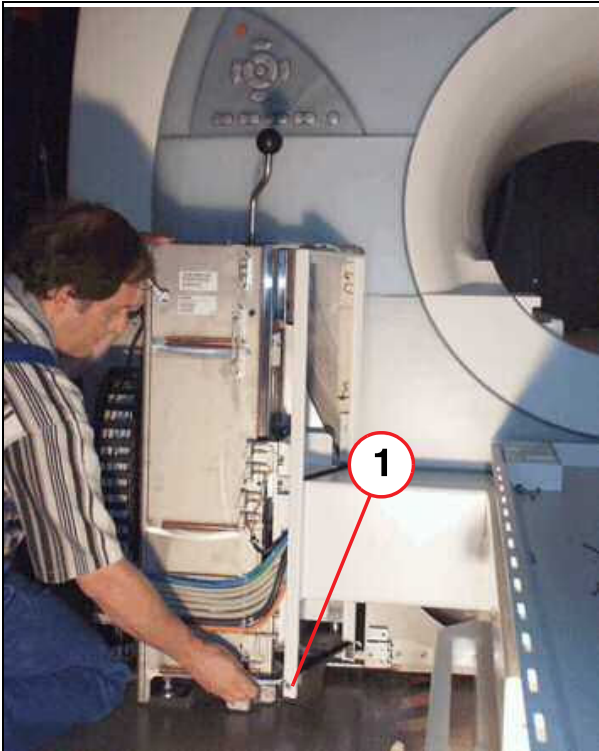

- Separate the bellows from the adhesive tape and remove the M6 Allen screw from the bellows guide.
- Use the crank and move the table down by hand.

Fig. 7: Remove the M6 Allen screws of the bellows guide



- Remove the M6 Allen screws of the bellows guide.

Fig. 8: Lower Allen screws of bellows guide



- Remove the lower M5 Allen screws.

Fig. 9: Remove hex bolt of bellows guide



- Remove hex bolt of bellows guide.
- Remove the bellows guide after removing all screws.

3.1.4.3 Support Beam

Fig. 10: Remove M8 Allen screw



- Remove the four M8 Allen screws, two above and two below.

Fig. 11: Remove two Allen screws inside the support beam



- Remove two Allen screws inside the support beam.
- Move the table up by hand, use the crank and mount the table service device.

- Mount the table transport device (service aid) (→ Fig. 12 Page 22).
 - Rollers must be evenly positioned/must be level on the floor so that you can return to the guide attachment position (→ Fig. 13 Page 22).



CAUTION

Possible damage of the spindle nut.


- » After the four transport rollers reached the floor and the load is removed from the spindle nut, do not continue to lower the table with the manual crank

- Lower the table (using the crank) (→ Fig. 14 Page 22) until the patient table is resting on the table service device. After the four transport rollers reached the floor and

the load is removed from the spindle nut, do not continue to lower the table with the crank.

- It takes 5 to 10 turns of the crank to remove the load from the spindle nut located behind the lifting carriage plate!

Fig. 12: Table transport device installation

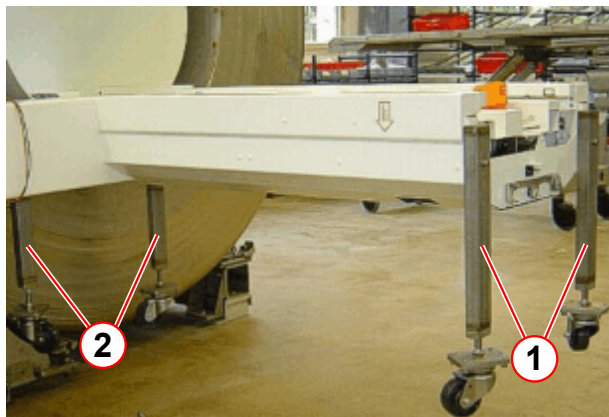

- (1) Transport roller 1
- (2) Transport roller 2

Fig. 13: Transport roller



Fig. 14: Lower the table with the crank until the load is removed from the spindle nut.

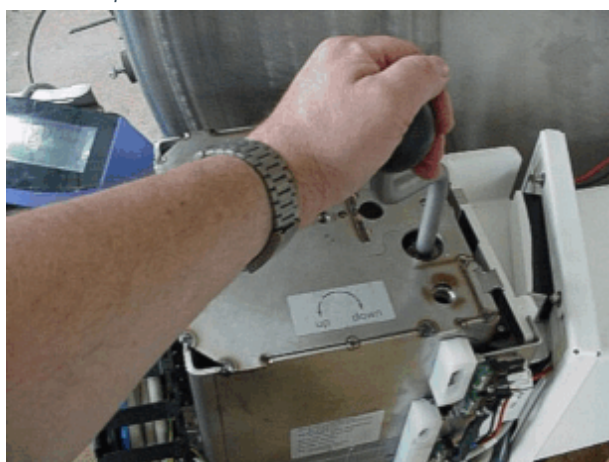


Fig. 15: Energy chain for lifting column




- Dismount the energy chain.

Fig. 16: Remove the M6 Allen screw



- Remove the M6 Allen screw of the spindle nut driver.


Fig. 17: M10 Allen screw, inside the support beam

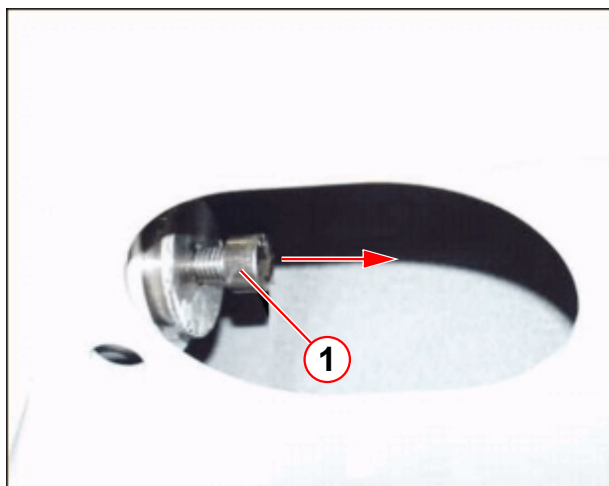

(1) M10 Allen screw

- Pull back the M10 Allen screw (inside the support beam) to release the spindle nut from the spindle nut driver. The table with the lifting column is now free mechanically!


Fig. 18: M8 Allen screw of support beam

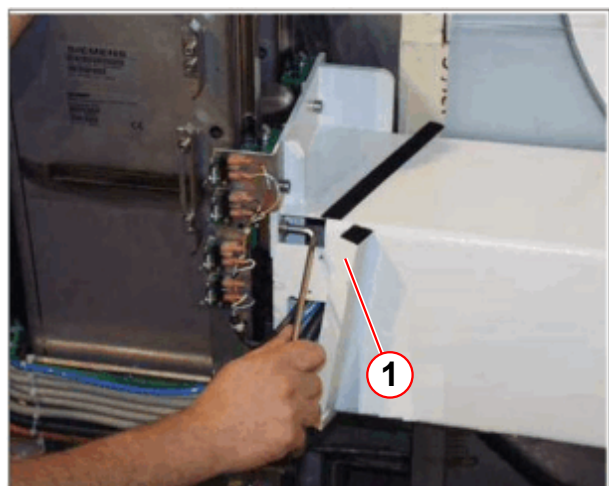

(1) Bellow support

- Remove the bellows support and remove the remaining M8 Allen screws.

Fig. 19: Remove the rest of all (16 in total) other M8 Allen screws

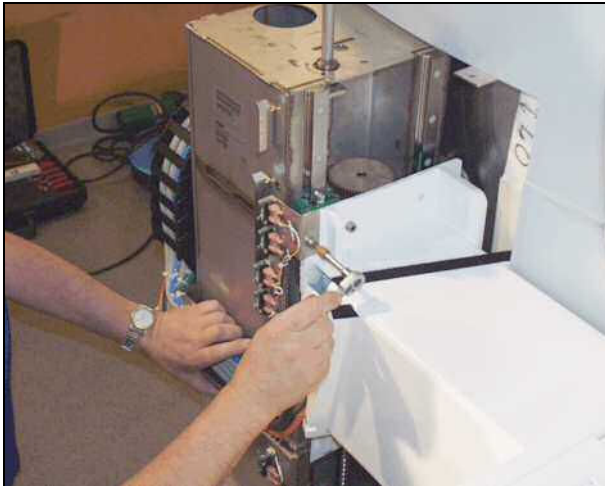

- Remove all other M8 Allen screws and carefully remove the switch panel.
- To access the guide rails, carefully push aside the table.


3.1.4.4 Guide rails

Fig. 20: Plastic covers from the guide rails



- Remove closing plugs if necessary (only used until 06/2008) by using an electric drill with a 5 mm bit and bore up the plastic closing plugs.

Fig. 21: Clinch rings

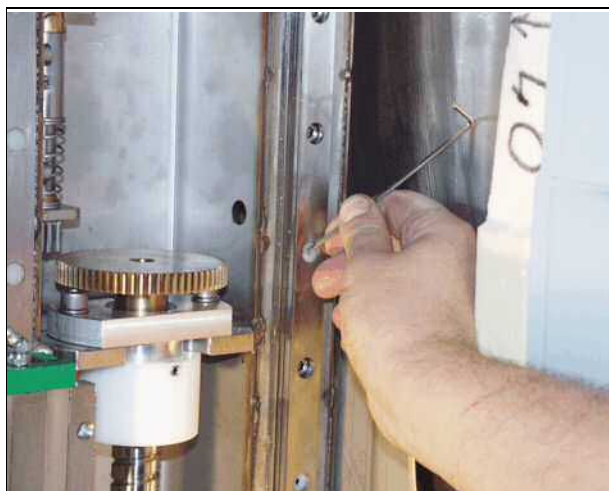


- Remove plastic clinch rings if necessary using a small screwdriver.

- Remove all M6 Allen screws and remove guide rail.

Fig. 22: Allen screws of guide rails



Fig. 23: Guide rail, removed

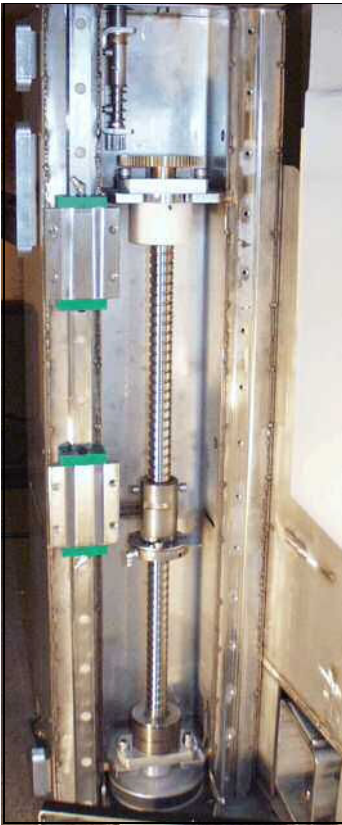

- Repeat the disassembly procedure on the other guide rail.

3.1.5 Assembly

3.1.5.1 Guide rails




Take into consideration the right and left side of the guide rail when mounting the new ones (ground edges of guide cassettes facing inward (→ Fig. 25 Page 28)).

- The two bearings remain at the bottom of the guide rail.

Fig. 24: Tighten the replaced guide rail from the middle up



Fig. 25: Ground edges of guide cassettes

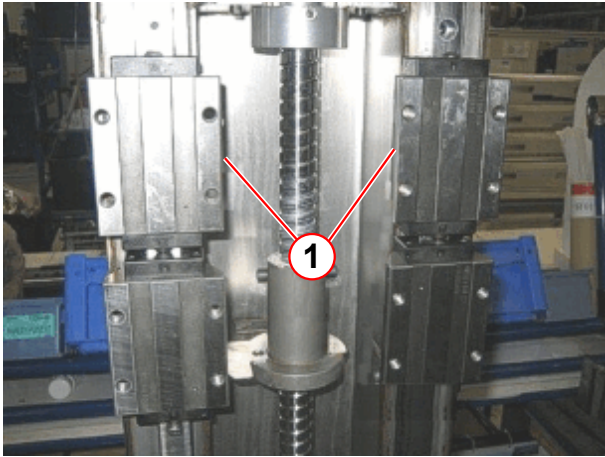


(1) Ground contact surface

- Reinsert Allen screws (apply Loctite 221) into counterbores from the upper part of the guide rail and tighten them. The bearings remain at the bottom of the guide rail.

- Move the bearings up. Insert one Allen screw as far as the bearings can remain on it.

Fig. 26: Move bearings up, insert screw, to bear on the screw

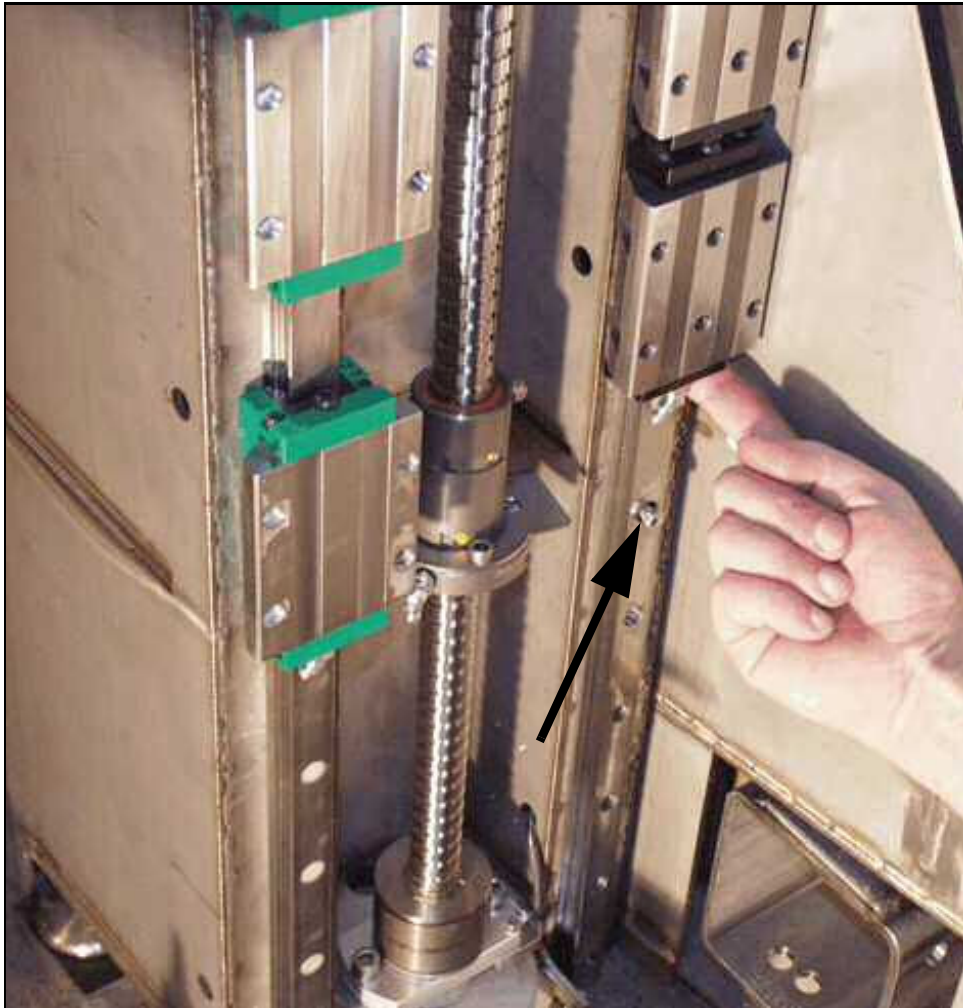

- Reinsert the rest of the Allen screws (apply Loctite 221) into counterbores from the lower part of the guide rail and tighten them.



The bearings must not drop down with force, as this would break the plastic end covers. Always be careful when moving the bearings.

- Repeat the procedure with the other guide rail.
- Before beginning with the reassembly, make sure the bearings are mounted correctly, the guide rails are clean, and no objects are present.


3.1.5.2 Support beam

Fig. 27: Move table back



- Carefully move the table back to the lifting column.
- Move the bearings up and install them on the support beam with one screw on each corner.
- The bearings must be seated accurately to the support beam.
- There are 4 fixing holes each for correct mounting.

Fig. 28: M8 Allen screw of support beam

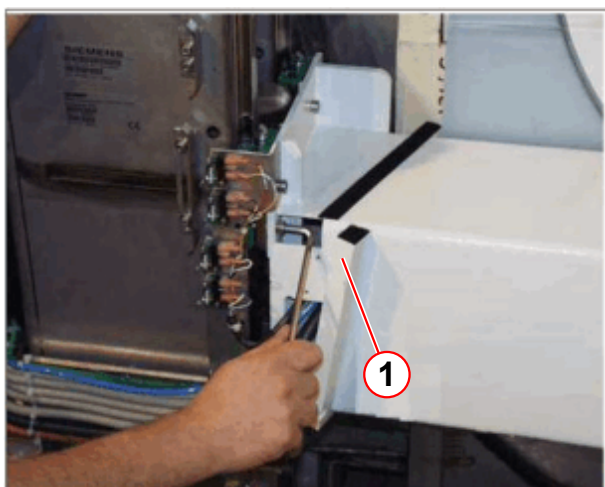

- Reassemble the switch bracket and insert all other Allen screws (apply Loctite 221) and tighten them.

(1) Bellow support

Once the support beam has been installed, the spindle nut in the column must be threaded back through the nut driver.


- Turn the spindle nut by hand so that it moves up underneath the nut driver and hold it in a diagonal/slanting position (approx. 45°). Collision switch S30 visible from underneath (→ Fig. 29 Page 31), (→ Fig. 30 Page 31).
- Use the crank to continue moving the spindle nut up so that it is inserted far enough into the nut driver.
- Turn the spindle nut in the nut driver back to the correct position (→ Fig. 31 Page 31). This ensures that the semicircular switching plate is in the right position to be able to actuate the cam switch of collision switch S30.

- Hold the spindle nut in the correct position and push the nut driver back on the M10 (→ Fig. 17 Page 24).

Fig. 29: Spindle nut




Fig. 30: Positioning

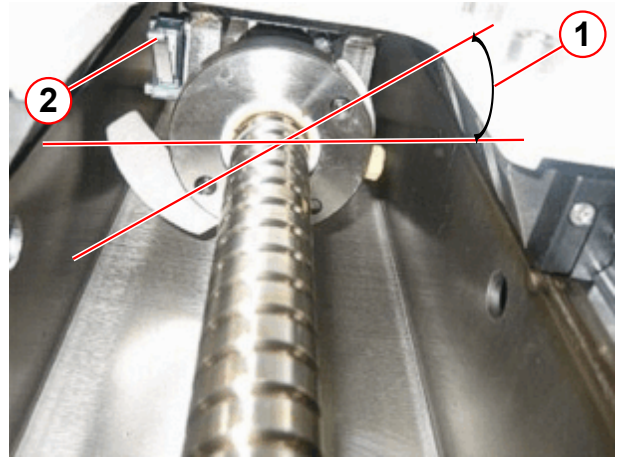


- (1) Approx. 45°
(2) Collision switch S30

Fig. 31: Correct position

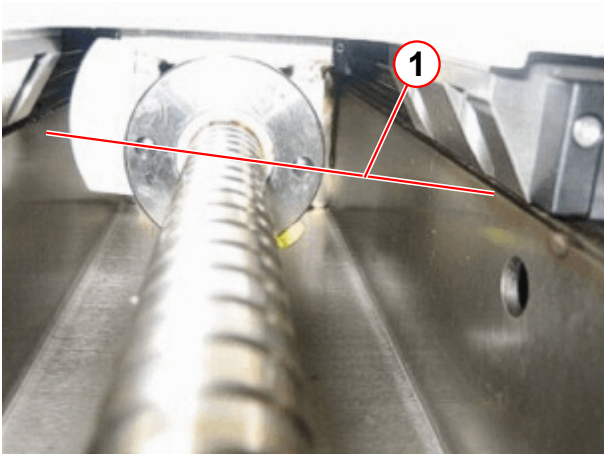

- (1) Correct position

Fig. 32: M10 Allen screw, inside the support beam

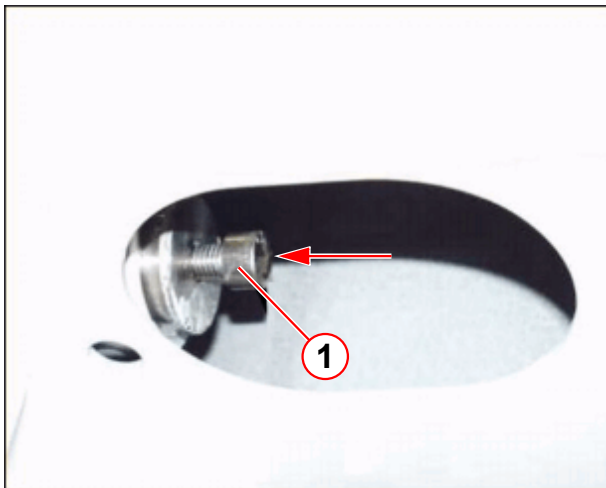

() M10 Allen screw

- Pull back the M10 Allen screw (inside the support beam) to connect the spindle nut to the spindle nut driver.

- Tighten the M6 Allen screw of the spindle nut driver.

Fig. 33: M6 Allen screw of spindle nut driver





Before applying the load again, the nut driver should display a small amount of end play and rotational play.

The table load is applied again by continuing to crank upwards.

This relieves the load on the transport device.

3.1.6 Concluding steps



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

Fig. 34: Energy chain for lifting column



- Reassemble the energy chain.
- Move the table up by hand, use the crank and remove the table service device.
- Install all covers in reverse order.
- Switch ON F22 (RF-room) in the ACC.

3.1.7 Adjustments

- n.a.

3.1.8 Final checks

- Switch on the patient table and check for proper function.
- Perform the safety test of the S30 spindle nut safety switch
- Ramp up the magnet and do a check of the magnetic field.

3.2 Lifting motor / toothed belt of vertical drive

The vertical lifting motor consists of magnetic objects. When replacing the motor, use caution around the magnetic field. Maintain a safe distance to the magnetic field during this procedure.



CAUTION

The vertical motor consists of magnetic objects.

Magnetic objects may cause injury.

- » Maintain a safe distance to the magnetic field.

3.2.1 Tools and time required

3.2.1.1 Time

Approx. 30 min.

3.2.1.2 Tools

Non-magnetic tool kit

3.2.2 Preliminary

- Switch OFF F22 (RF-room).

3.2.3 Disassembly

3.2.3.1 Covers



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

- Remove the service access cover and the lifting column.

The removal of the service access cover and the lifting column cover is described in the REP/system, refer to (→ System / M6-020.841.10).

3.2.3.2 Vertical motor



Drive the table to the lowest position before replacing the motor or the tooth belt.

Fig. 35: Remove covers and disconnect cables



- Remove the patient table electronics plate and disconnect the cables belonging to the vertical motor.

Fig. 36: Release toothed belt




- Release toothed belt, press on the motor to remove toothed belt, and replace belt, if necessary.



The vertical lifting motor consists of magnetic objects. When replacing the motor, use caution around the magnetic field. Maintain a safe distance to the magnetic field during this procedure.

Fig. 37: Release M12 nuts and remove motor

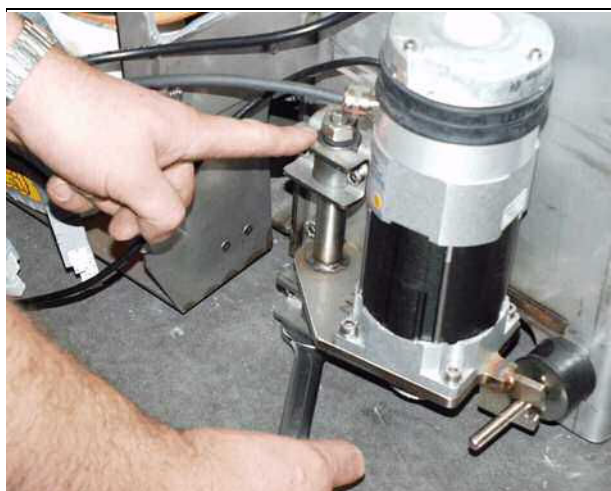

- Remove M12 nuts and turn the motor aside.
- Remove the four Allen screws and **remove the motor immediately from the magnetic field.**

3.2.4 Assembly



The vertical lifting motor consists of magnetic objects. When replacing the motor, use caution around the magnetic field. Maintain a safe distance to the magnetic field during this procedure.

- Replace motor and safety immediately with the four Allen screws. **When replacing the motor, use caution around the magnetic field.**
- Move the motor back into its original position.
- Reinsert the toothed belt.
- Stretch the toothed belt and tighten the M12 nuts.

3.2.5 Concluding steps

- Reconnect all cables.
- Switch ON F22 (RF-room).
- Install covers.

3.2.6 Adjustments

- n.a.

3.2.7 Final checks

- Switch on the patient table and verify proper vertical movement.

3.3 Spindle

3.3.1 General

The table is moved in vertical direction via a special bearing spindle. That spindle consists of a toothed wheel, upper spindle bearing, end stop, spindle nut, lower spindle bearing and lower toothed wheel.



There are two different methods described for replacing the spindle.

If the main bearing is defective it may be necessary to use **method I** with the **table transport device**.

If the vertical movement is still possible, however, the spindle should be replaced due to a damaged brake, use **method II** with a **massive wooden** beam. The wooden beam is part of the spare parts delivery.

The spindle can be replaced by one person.



CAUTION

Strong magnetic field!

Ferromagnetic spindle may cause injury.

» Ramp down the magnet




Two types of spindles are introduced. The newer type has an additional bearing on the spindle brake.

In some pictures, minor differences may occur.


New version of spindle bearing

Fig. 38: Vertical motor

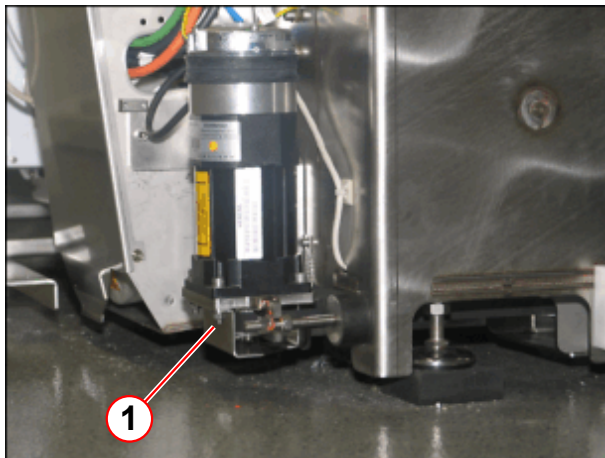


(1) Additional safety cover

Fig. 39: New version of brake bearing

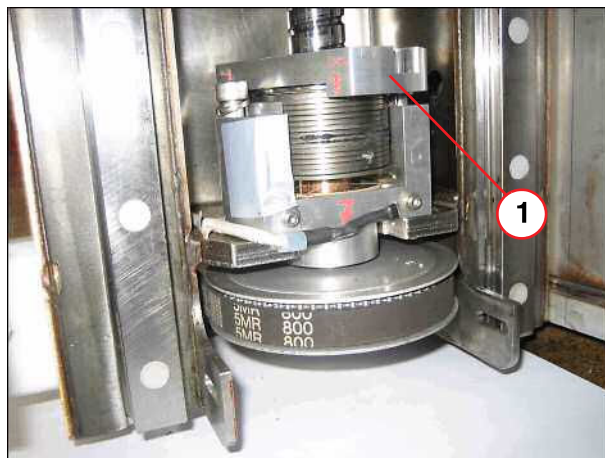


(1) Additional bearing

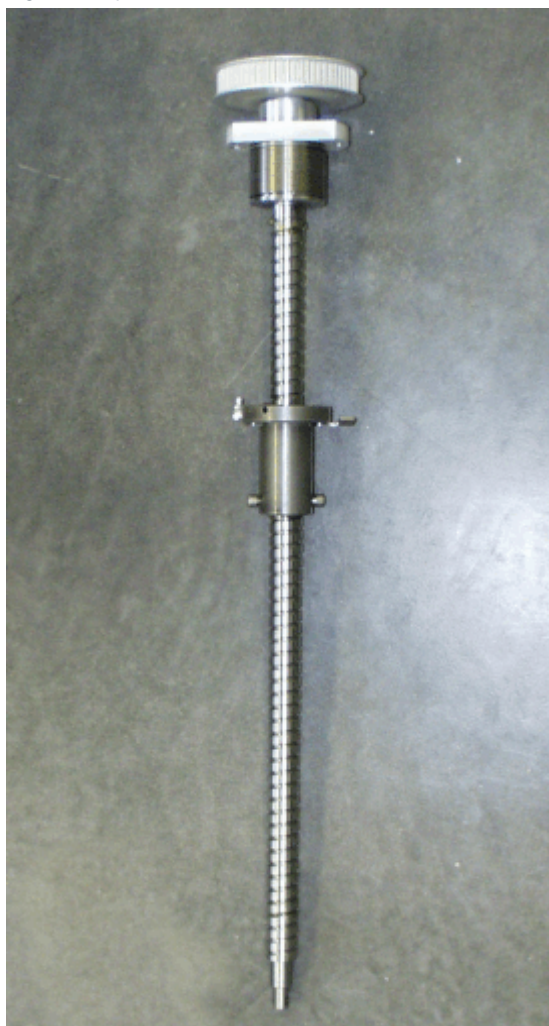Fig. 42: Spindle*Fig. 43: Spare material*

3.3.3 Principle of spindle (old version)

Fig. 44: Principle of spindle

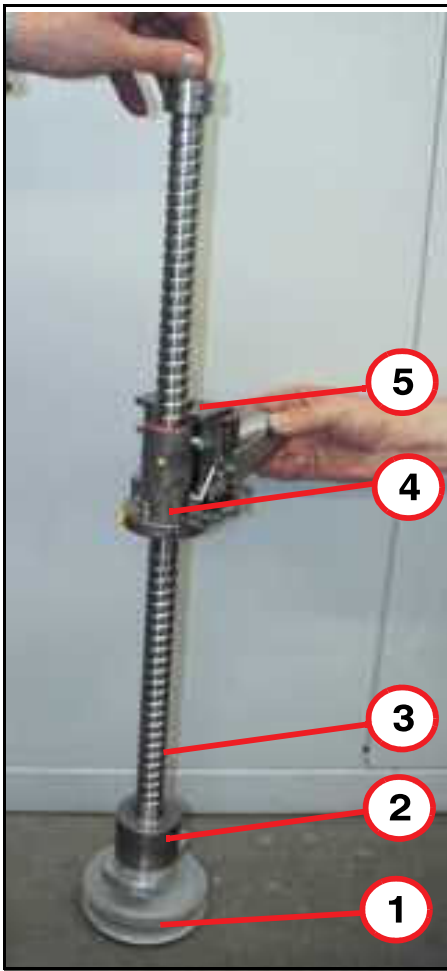

- (1) Toothed wheel, lower spindle
- (2) Mechanical brake
- (3) Spindle:
- (4) Spindle nut:
- (5) Spindle nut driver

The spindle consists of:

- Lower spindle toothed wheel, driven by a horizontal motor (1).
- Mechanical brake, active during down movement (2).
- Spindle nut with ball bearing circulation spindle (4).
- Spindle nut driver (5).

3.3.4 Tools and time required

3.3.4.1 Time

The estimated time does **not include** ramping the magnet.


Method I, with table transport device, approx. 120 min

Method II, with wooden beam, approx. 60 min

3.3.4.2 Tools

Table transport device (method I).

Fig. 45: Table transport device

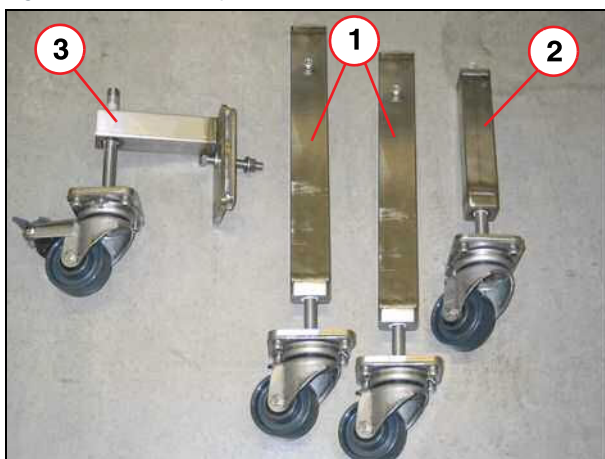

- (1) Transport roller 1
- (2) Transport roller 2
- (3) Transport roller 3

Massive wooden beam as part of the spare parts delivery for the spindle.

Non-magnetic tool kit.

Loctite 221.

Working gloves.

Ratchet and socket kit.

3.3.5 Preliminary

- Ramp down the magnet



Any use of grease or oil for the patient table spindle may reduce the brake force and will cause a vertical brake fault.

3.3.6 Disassembly method I with table transport device

3.3.6.1 Covers



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

- Remove the service access cover and the lifting column cover.

The removal of the service access cover and the lifting column cover is described in the REP/system, refer to (→ System / M6-020.841.10).

Fig. 46: Lower magnet, front cover



- Remove lower front cover of the magnet.

3.3.6.2 Bellows guide

Fig. 47: Remove the M6 Allen screws of the bellows guide



- Separate the bellows from the adhesive tape and remove the M6 Allen screw from the bellows guide.
- Remove the M6 Allen screws of the bellows guide.

Fig. 48: Lower Allen screws of bellows guide



- Remove the lower M5 Allen screws.

Fig. 49: Remove hex bolt of bellows guide



- Remove the M6 Hex bolt of the bellows guide.

Fig. 50: Lower bellows guide/toothed wheel




- Separate the bellows from the adhesive tape and remove the M6 Allen screw from the bellows guide.
- Remove the bellows guide after removing all screws.
- Remove lower bellows bracket

3.3.6.3 Support beam

- Move the table up by hand, use the crank.
- Mount the table service device and lower the table until the lifting column rests on the table service device.
- Switch OFF F22 (RF- room) in the ACC.

Fig. 51: Table transport device installation

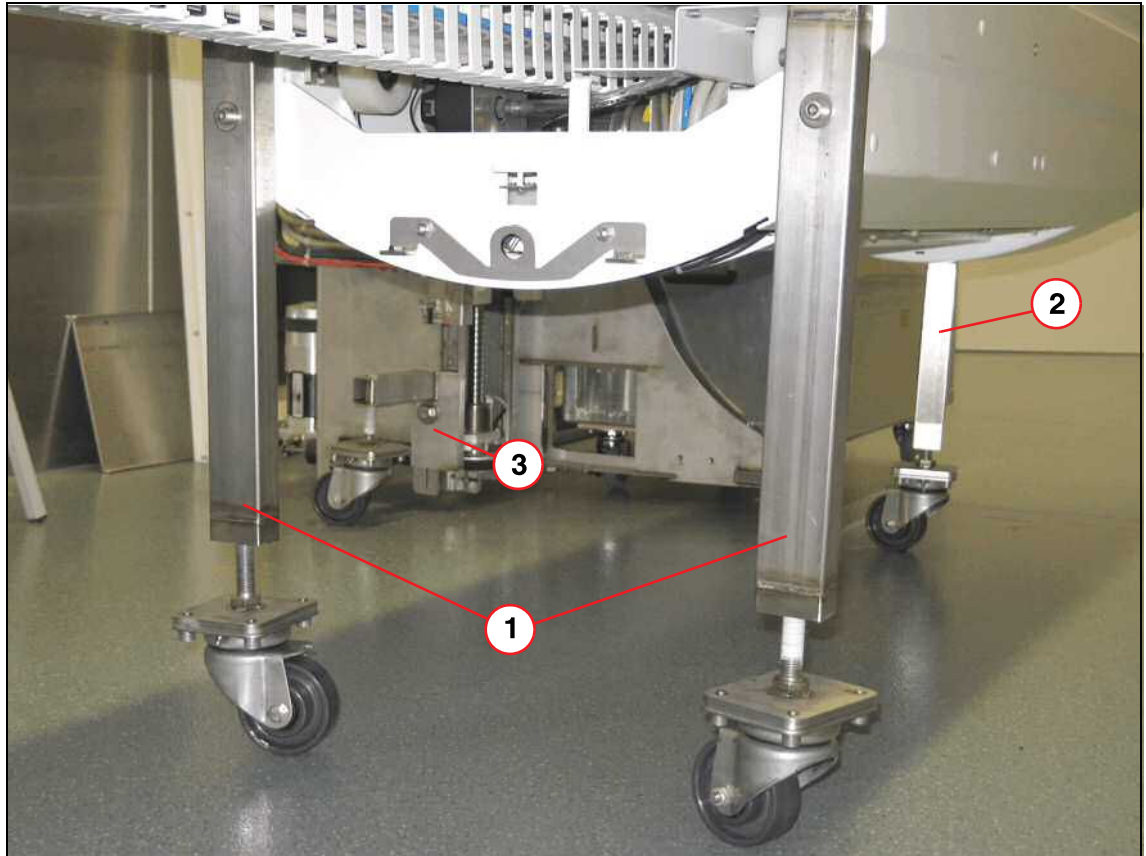

- (1) Transport roller 1
- (2) Transport roller 2
- (3) Transport roller 3

Fig. 52: Remove M8 Allen screw



- Remove the four M8 Allen screws, two above and two below.

Fig. 53: Remove the M6 Allen screw




- Remove the M6 Allen screw of the spindle nut driver.

Fig. 54: M6 Allen screw, inside the support beam



- Pull back the M6 Allen screw (inside the support beam) to release the spindle nut from the spindle nut driver. The table with the lifting column is now mechanically separated from the spindle.

Fig. 55: M8 Allen screw of support beam

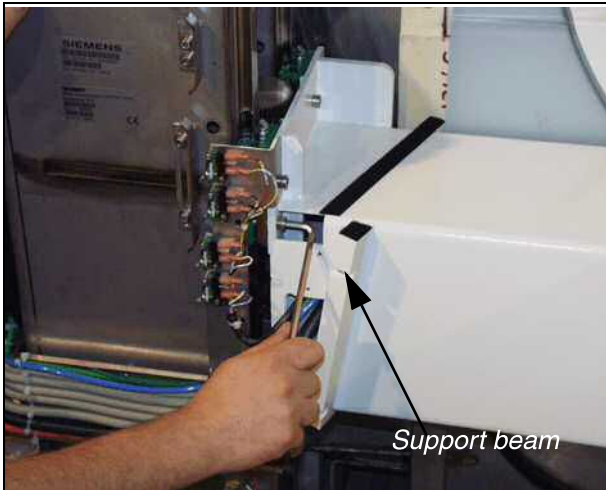

- Remove the bellows support and remove the remaining M8 Allen screws.

Fig. 56: Remove the rest of all (16 in total) other M8 Allen screws



- Remove all other M8 Allen screws and carefully remove the switch panel. Make sure the lifting column rests safely on the table service device.
- To access the spindle, carefully push the table to the side.

3.3.6.4 Toothed belt

Fig. 57: Release toothed belt



- Release toothed belt, press on the motor to remove toothed belt, and replace belt, if necessary.



If no motor brake is active, the table can **move down** slowly.

3.3.6.5 Cog wheel

Fig. 58: Upper toothed wheel



(1) Grub screw of upper toothed wheel

- Release the grub screw and remove the cog wheel.

3.3.6.6 End stop

Fig. 59: Upper mechanical end stop



(1) Upper mechanical end stop removed

- Remove the mechanical end stop.

3.3.6.7 Guide bearing

Fig. 60: Spindle bearing



(1) Upper spindle bearing

- Remove both M8 Allen screws from the upper guide bearing.
- If the upper guide bearing cannot be removed slightly, use the longer screws delivered with the spare part.

3.3.6.8 Spindle

Fig. 61: Lower spindle component



- (1) Temperature sensor
(2) Lower spindle, bearing block

- Remove the M3 Allen screw for the temperature sensor and both M8 Allen screws for the lower spindle bearing.

- Remove the spindle and replace it. Make sure the thread from the temperature sensor is positioned properly.



If there are washers bonded to the bottom of the spindle bearing, don't remove them.

Fig. 62: Remove spindle



3.3.7 Disassembly method II, with a massive wood beam

3.3.7.1 Covers



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.



For the vertical movement during exchange procedure II, the patient table can be driven with the normal motor function. Use caution when moving the patient table.

- Remove the service access cover and the lifting column cover.

The removal of the service access cover and the lifting column cover is described in the REP/system, refer to (→ System / M6-020.841.10).

Fig. 63: Lower magnet, front cover



- Remove lower front cover of the magnet.

3.3.7.2 Bellows guide

Fig. 64: Remove the M6 Allen screws of the bellows guide



- Separate the bellows from the adhesive tape and remove the M6 Allen screw from the bellows guide.
- Remove the M6 Allen screws of the bellows guide.

Fig. 65: Lower Allen screws of bellows guide



- Remove the lower M5 Allen screws.

Fig. 66: Remove hex bolt of bellows guide



- Remove the M6 Hex bolt of the bellows guide.

Fig. 67: Lower bellows guide/toothed wheel



- Separate the bellows from the adhesive tape and remove the M6 Allen screw from the bellows guide.
- Remove the bellows guide after removing all screws; if necessary, lower the table.
- Remove lower bellows bracket.

3.3.7.3 Cog wheel

Fig. 68: Upper toothed wheel



(1) Grub screw of upper toothed wheel

- Drive the table down until the upper spindle bearing and the toothed wheel are accessible.
- Release the grub screw and remove the cog wheel.

3.3.7.4 End stop

Fig. 69: Upper mechanical end stop



(1) Upper mechanical end stop removed

- Remove the end stop.

3.3.7.5 Guide bearing

Fig. 70: Spindle bearing




(1) Upper spindle bearing

- Remove both M8 Allen screws from the upper guide bearing.
- If the upper guide bearing cannot be removed slightly, use the longer screws delivered with the spare part.

3.3.7.6 Spindle


Fig. 71: Principle of spindle

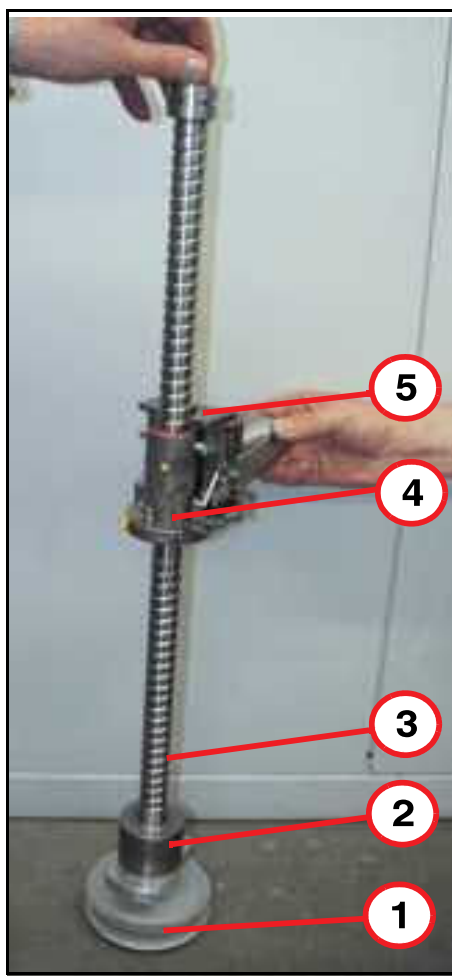

- (1) Toothed wheel, lower spindle
- (2) Mechanical brake
- (3) Spindle:
- (4) Spindle nut:
- (5) Spindle nut driver

- Move the table up by pressing the table **UP** button. The lower edge of the table (support beam) must reach a height of approx. 600 mm from the floor.

Fig. 72: Support frame

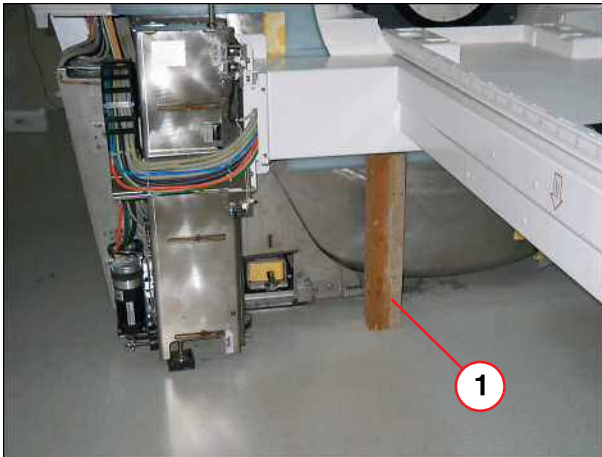

(1) Wooden beam under support frame

- Use a massive wooden beam and place it under the support frame. Wooden beam delivered with the spindle spare part. (Alternative use a frame with a height of 550 mm to carry a load of 220 kg).
- Lower the table by pressing the table **DOWN** button until the lifting column rests on the massive wooden beam or frame. The error (MRI_PXX) ID **63** will show up.
- Switch OFF F22 (RF-room) in the ACC.

Fig. 73: Release toothed belt



- Release the nut and push the motor to remove the toothed belt.



If no motor brake is active, the table can **move down** slowly.

Fig. 74: Remove the M6 Allen screw




- Remove the M6 Allen screw of the spindle nut driver.

Fig. 75: M6 Allen screw, inside the support beam



- Pull back the M6 Allen screw (inside the support beam) to release the spindle nut from the spindle nut driver. The table with the lifting column is now free mechanically!

Fig. 76: Function of spindle nut

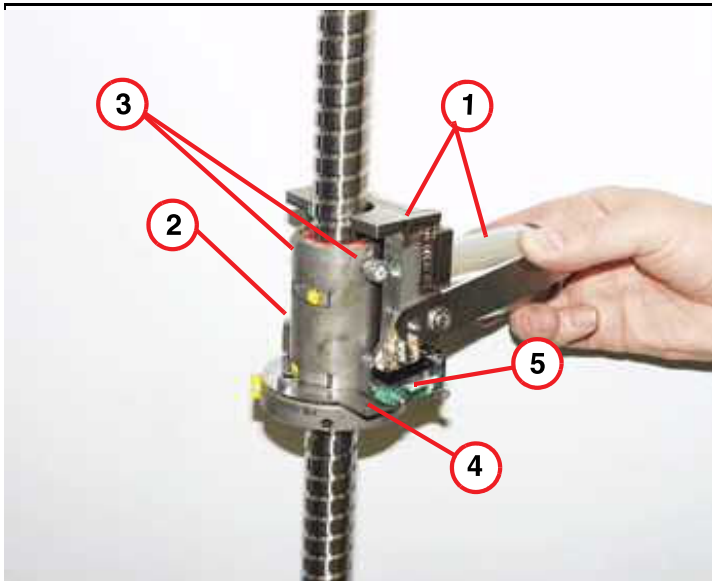


- (1) Spindle nut driver
- (2) Spindle nut:
- (3) Allen screws
- (4) Switchplate
- (5) Switch S30

- The principle behind disassembly the spindle: Is to turn the spindle nut without de-mounting the support frame so that the Allen screws are no longer in contact with the spindle nut driver.

As a result, the spindle nut has to be turned.


Fig. 77: Spindle nut prepared for spindle replacement

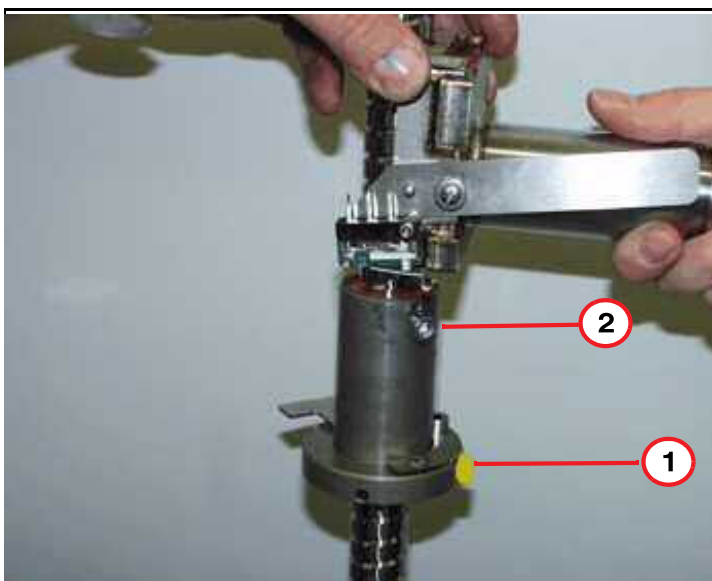

- (1) Turned lubrication nipple
- (2) Turned Allen screws of spindle nut

- To release the spindle nut from the spindle nut driver: Turn the spindle nut until the lubrication nipple looks towards the support beam, while the spindle nut remains in this position.
- The spindle nut can be driven down by holding the spindle nut and turning the lower toothed wheel by hand.

Fig. 78: Principle of spindle nut replacement

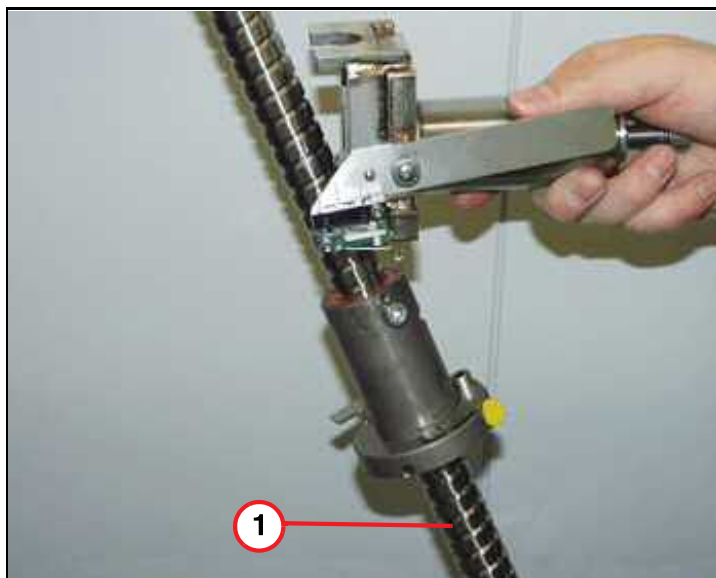

(1) Spindle tilt to weave it out

- In case the spindle nut gets stuck due to the spindle nut driver, tilt the top of spindle. In this way, the spindle has more room to pass the spindle nut while driving down.

Fig. 79: Lower spindle component



(1) Temperature sensor
(2) Lower spindle, bearing block

- Remove the M3 Allen screw for the temperature sensor and both M8 Allen screws for the lower spindle bearing.



If there are washers bonded to the bottom of the spindle bearing, don't remove them.

- Remove the spindle and replace it. Make sure the thread from the temperature sensor is positioned properly.

Fig. 80: Remove spindle



3.3.8 Assembly method I

3.3.8.1 Spindle

- Mount the new spindle with both M8 Allen screws and safety the lower bearing, re-connect the temperature sensor.
- Position the spindle nut by turning the cog wheel by hand (until the spindle nut reaches the same level as the support beam).
- Tighten the M8 Allen screws of the lower and upper spindle bearing . Make sure the thread from the temperature sensor is positioned properly.



If no motor brake is active, the table can **move down** slowly.

- Insert the upper cog wheel and tighten the grub screw.
- Insert the end stop.

3.3.8.2 Support beam

Fig. 81: Move table back

- Carefully move the table back to the lifting column.
- Move up the bearings and install them to the lifting column with one screw at each corner.

Fig. 82: Safety bearings with Allen screws

- Reassemble the switch bracket, insert all Allen screws and tighten them.

Fig. 83: M6 Allen screw, inside the support beam



- Turn the spindle nut into position and push the M6 Allen screw (inside the support beam) back to contact the spindle nut driver.

- Insert the M6 Allen screw (apply Loctite 221) and safety the spindle nut driver.

Fig. 84: M6 Allen screw of spindle nut driver




Make sure the M6 Allen screw of the spindle nut driver **cannot** be pulled back.

3.3.9 Assembly method II

3.3.9.1 Spindle


Fig. 85: Spindle nut prepared for spindle replacement

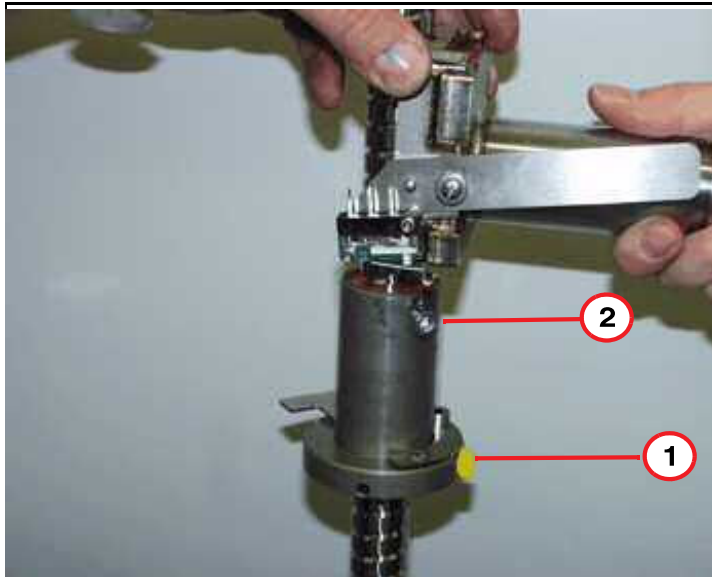

- (1) Turned lubrication nipple
- (2) Turned Allen screws of spindle nut

- Insert the new spindle and safety the lower spindle bearing with both M8 Allen screws. Reconnect the temperature sensor. Make sure the thread from the temperature sensor is positioned properly.
- Drive the spindle nut up by holding it and turning the lower toothed wheel by hand. Until the spindle nut reaches the same level as the spindle nut driver.
- Turn the spindle nut until the lubrication nipple looks toward the support beam, while the spindle nut remains in this position.
- Drive the spindle nut up into the spindle nut driver: By holding the spindle nut in this position and turning the lower toothed wheel by hand until it stops.

Fig. 86: Function of spindle nut

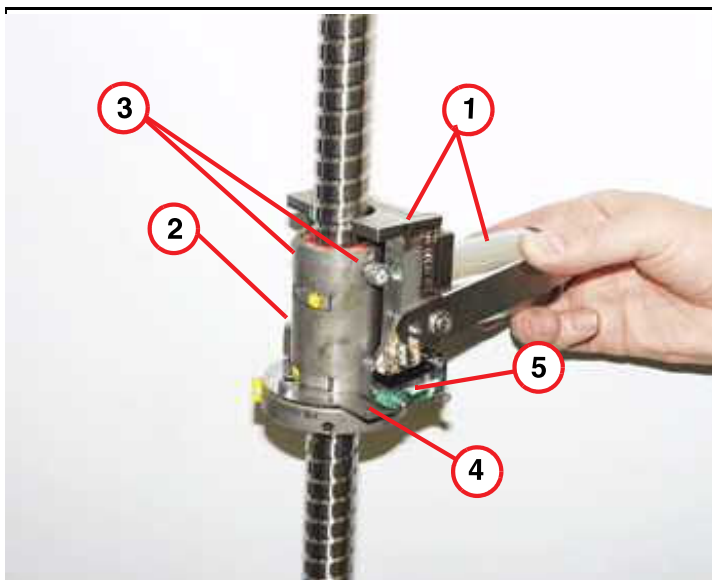

- (1) Spindle nut driver
- (2) Spindle nut:
- (3) Allen screws
- (4) Switchplate
- (5) Switch S30

- Turn the spindle nut and lower the spindle nut until both Allen screws are positioned in the spindle nut driver. The switch bracket of S30 is in the correct position.

Fig. 87: M6 Allen screw, inside the support beam



- Turn the spindle nut into position and push the M6 Allen screw (inside the support beam) back to contact the spindle nut driver.

- Insert the M6 Allen screw (apply Loctite 221) and safety the spindle nut driver.

Fig. 88: M6 Allen screw of spindle nut driver



Make sure the M6 Allen screw of the spindle nut driver **cannot** be pulled back.

- Switch on F22 (RF- room) in the ACC.
- Drive the table up, by pressing the table **UP** button and remove the wooden frame, lower the table until the upper spindle part is accessible.



If no motor brake is active, the table can **move down** slowly.

- Remount the upper spindle bearing and tighten the M8 Allen screws of the lower and upper spindle bearing.
- Insert the upper cog wheel and tighten the grub screw.
- Insert the end stop.

3.3.10 Concluding steps

- Install the toothed belt of the vertical drive and tighten it. If necessary, replace belt.
- Install the lower bellows bracket.
- Install the bellows guide.
- If necessary switch ON F22 (RF-room) ptab and verify proper vertical movement.
- Drive the table up and down (approx. 20 times) until the spindle is getting hot and the error message **(MRI_PXX) 66** shows up. Wait until the spindle is cooled down.
- Install the covers.

3.3.11 Adjustments

- n.a.

3.3.12 Final checks

- Verify proper vertical movement.

3.4 Brake disk

3.4.1 General

The brake disk is responsible for the brake torque during downward movement of the patient table. The brake disk can be replaced (in accordance with the serial numbers) with the magnet at field as long as anti-magnetic tools are used.



CAUTION

Strong magnetic field.

Magnetic objects may cause injury.

- » **Maintain a safe distance to the magnetic field.**




Two types of spindles are introduced. The newer type has an additional bearing on the spindle brake.

In some pictures, minor differences may occur.


New version of spindle bearing

Fig. 89: Vertical motor

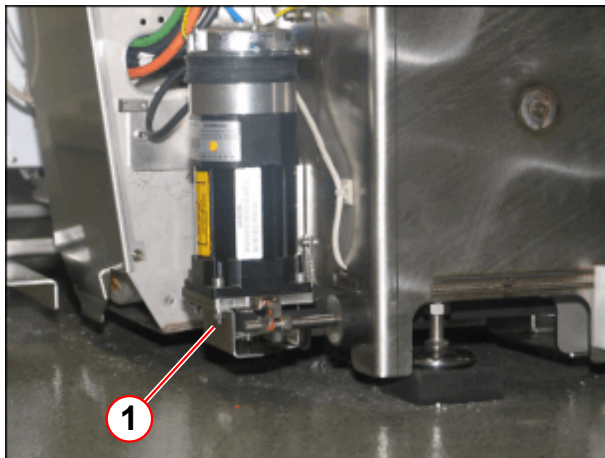

(1) Additional safety cover


Fig. 90: New version of brake bearing

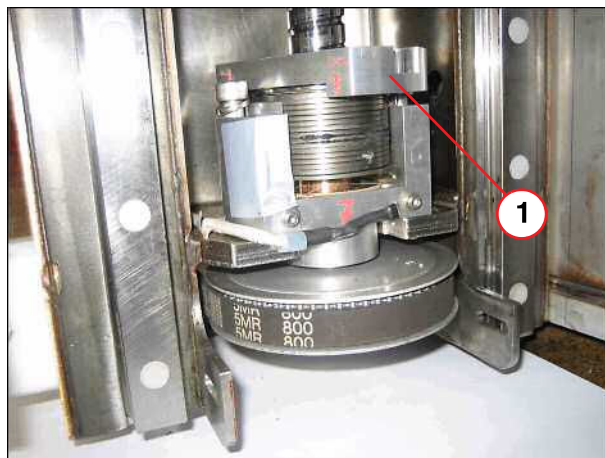

(1) Additional bearing

Old version of spindle bearing

Fig. 91: Toothed disk

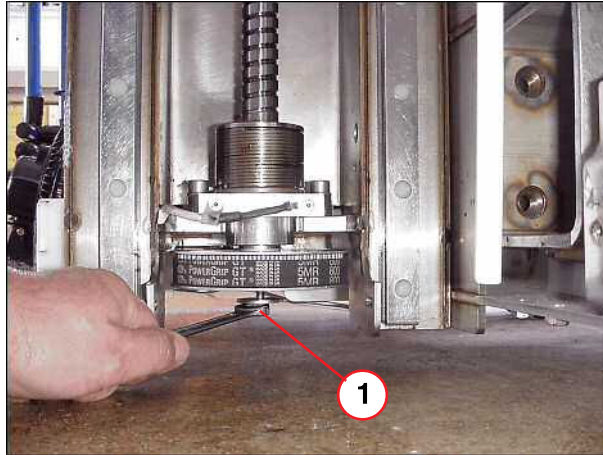

(1) Loose mounting screw of toothed wheel 10 mm to 12 mm

3.4.2 Tools and time required

3.4.2.1 Time

Approx. 60 min

3.4.2.2 Tools

Non-magnetic tool kit.

Wooden beam (delivered with the spare part disk).

3.4.3 Preliminary

- Only for basic tables with serial numbers > 1529, exchangeable tables with serial number > 1094, and all spindles replaced after June 2005
- A worn brake disk can be indicated by error 77 after the number of lifts exceeds 100000 movements.
- Switch OFF F22 (RF-room) in the ACC.

3.4.4 Disassembly

3.4.4.1 Covers



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

- Remove the service access cover and the lifting column cover.

The removal of the service access cover and the lifting column cover is described in the REP/system, refer to (→ System / M6-020.841.10).

Fig. 92: Lower magnet, front cover



- Remove lower front cover of the magnet.

3.4.4.2 Bellows guide

Fig. 93: Remove the M6 Allen screws of the bellows guide



- Separate the bellows from the adhesive tape and remove the M6 Allen screw from the bellows guide.
- Remove the M6 Allen screws of the bellows guide.

Fig. 94: Lower Allen screws of bellows guide



- Remove the lower M5 Allen screws.

Fig. 95: Remove hex bolt of bellows guide



- Remove the M6 Hex bolt of the bellows guide.

Fig. 96: Lower bellows guide/toothed wheel



- Separate the bellows from the adhesive tape and remove the M6 Allen screw from the bellows guide.
- Remove the bellows guide after removing all screws.
- Remove lower bellows bracket.

3.4.4.3 Brake disk




Any use of grease or oil for the patient table spindle may reduce the brake force and will cause a vertical brake fault.



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.


Fig. 97: Vertical motor

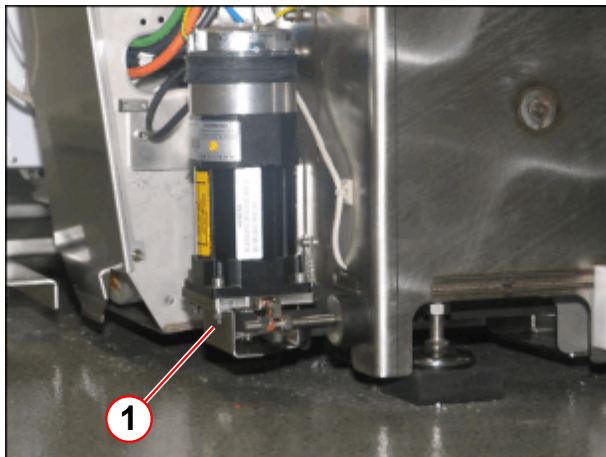

(1) Additional safety cover

- On newer system, remove the safety cover.


Fig. 98: Toothed disk

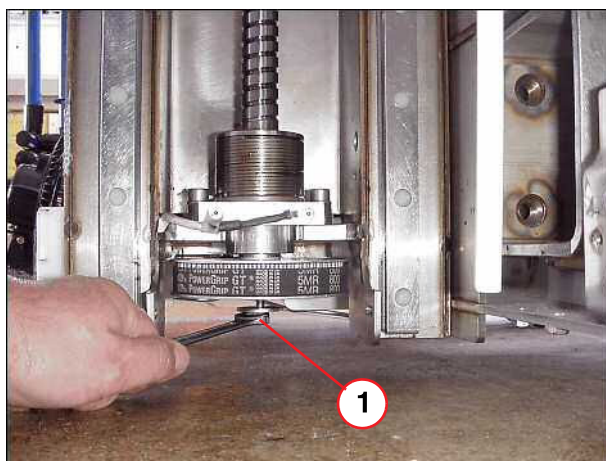

(1) Loose mounting screw of toothed wheel 10 mm to 12 mm

- Release mounting screw of toothed wheel approx. 10 mm- 12 mm (do not remove the screw).

Fig. 99: Support frame

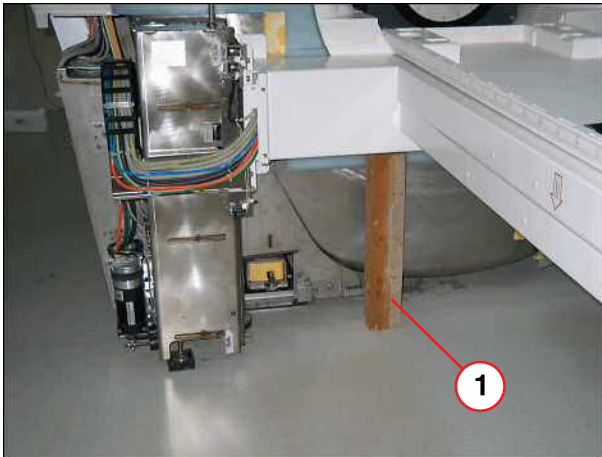

(1) Wooden beam under support frame

- Place the wooden beam under the support frame of the patient table.



CAUTION


Possible damage of the spindle nut.

- » After the wooden support beam is supporting the frame and the load is removed from the spindle nut, do not continue to lower the table with the manual crank.

Fig. 100: Crank for vertical movement



Fig. 101: Brake disk with clearance

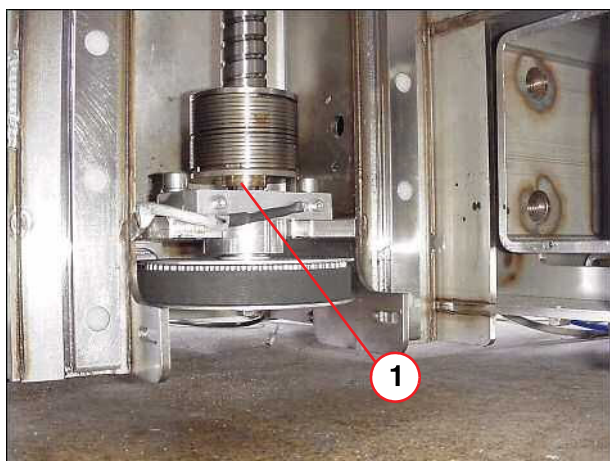

(1) Clearance between brake disk and mounting bracket approx. 10 mm - 12mm.

Fig. 102: Removing the brake disc



(1) Flat-bladed screw driver

- Crank the patient table down until the complete spring clutch is raised and a clearance of 10 mm- 12 mm is visible above the brake disk.

The clearance is limited by the released screw on the button of the toothed wheel.

!If the clearance is not visible, it might be necessary to lift up the spindle manually by using a big screwdriver!

After the wooden support beam is supporting the frame and load is removed from the spindle nut, do not continue to lower the table with the manual crank.

- Remove the old brake disk using a small flat bladed screwdriver.

Clean the clearance after removing the old brake disk.

Fig. 103: Lower surface: copper colored



Fig. 104: Upper surface: grooved



- A worn brake disk can be indicated by Error 77.

The color of the sliding surface will be primarily bronze or copper.

Clean the clearance after removing the old brake disk.



Make sure that the grooved side of the brake disc faces up.

Fig. 105: Brake disks



- The spare part disk halves are positioned below the clutch section.

Clean the clearance after removing the old brake disk.


- Using a small flat bladed screwdriver, first insert the half with the hole and close it in place with the securing pin.

Fig. 106: Replace new brake disk




(1) Flat blade screw driver

Fig. 107: Removing the brake disk (newer version)

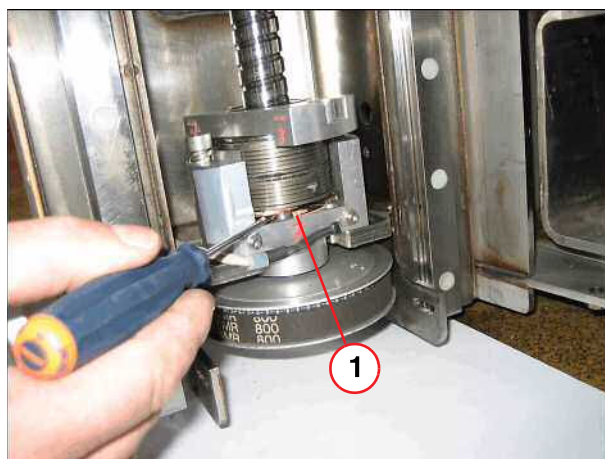

(1) Flat blade screw driver

Fig. 108: Brake disk

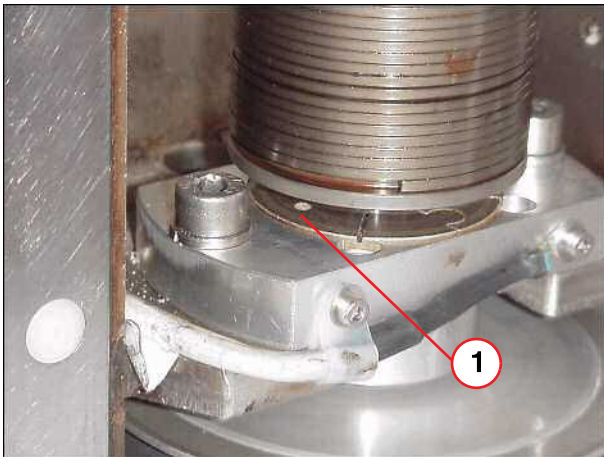

(1) One half of brake disk placed into the securing pin

- Both disk halves lie flush. The half with the hole is closed in place with the securing pin.

Both halves, one fitted into the other, have to lie flush.


Fig. 109: Crank for vertical movement



- Use the hand crank to move the table up until the spring clutch is in contact with the brake disk.
- Remove the wooden beam under the support frame.


There should be no clearance between the brake disk and the spring clutch.

Fig. 110: Brake disk in contact with spring clutch

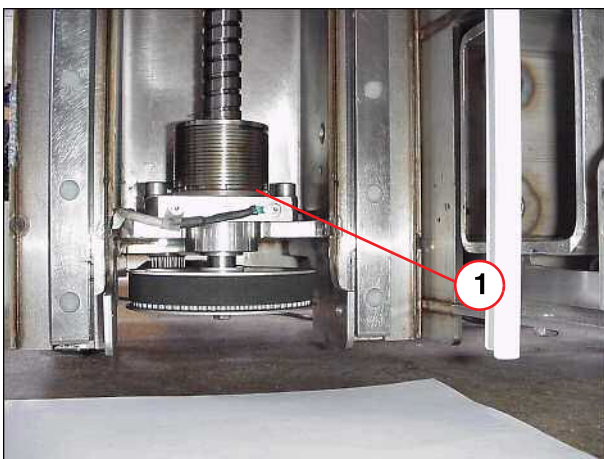


(1) Brake disk seated on the spring clutch

Fig. 111: Mounting screw of toothed disk

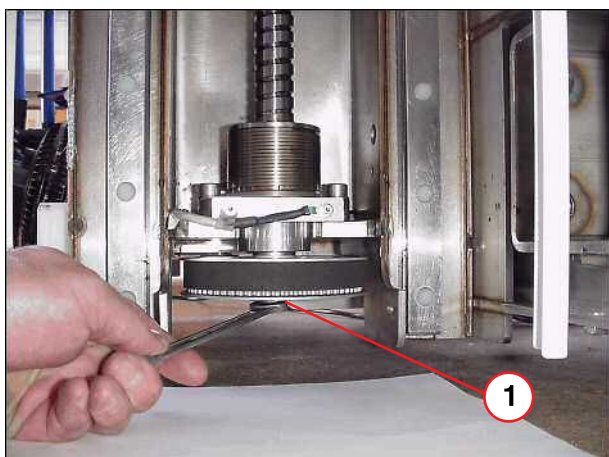

(1) Apply Loctite and tighten mounting screw

Fig. 112: Mounting screw of toothed disk (new version)

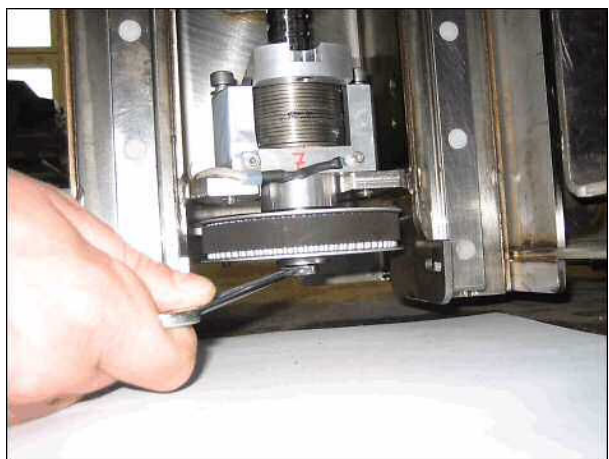

- Apply Loctite 221 and tighten the mounting screw of the toothed disk.
- If present, remount safety bracket.

3.4.5 Assembly

- Install the lower bellows bracket.
- Install the bellows guide.

3.4.6 Concluding steps

- Install the covers.

3.4.7 Final checks

- Switch ON F22 (RF-room) for the patient table and verify proper vertical movement.
- Initiate up and down movement several times. The spindle brake needs to be applied several times before full stopping power is attained.

3.5 Tabletop

3.5.1 General

When replacing the tabletop plate, we recommend a second person be used to help.

A minimum distance between the magnet and the wall is required for replacement of the tabletop plate. Depending on the system type, the minimum distances are:

- Avanto **1450 mm**
- Espree **1650 mm**
- Trio a Tim System **850 mm**
- Symphony a Tim System **1300 mm**
- Verio **1300 mm**

3.5.2 Tools and time required

3.5.2.1 Time

Approx. 60 min

3.5.2.2 Tools

Non-magnetic tool kit

3.5.3 Preliminary

- Switch OFF F22 (RF-room) in the ACC



According to the system/room specifications the movement range of the horizontal movement can be reduced.

With the whole body version of the horizontal drive, the movement range of the horizontal movement may be reduced. If that is the case, remove the end stops.
(→ Fig. 289 Page 173)

Before disassembly the patient table covers on a system with a partial drive unit: The toothed belt must be released to push the table top to the service end of the magnet. Refer to (→ Prerequisite for cover removal / partial drive / Page 98)

3.5.4 Disassembly

3.5.4.1 Covers



Use caution when moving the tabletop plate > 30 kg. A second person is necessary.



Be careful to avoid damaging any cables when removing the covers of the patient table. Use caution when moving the patient table.

Basic table

Fig. 113: Coil plug cover (foot end)



- Remove the cover of the coil plug.


Fig. 114: Remove Allen screws



- Remove the lower Allen screws of the cover from the coil plugs.

Removable table

Fig. 115: RF-Cover RTT

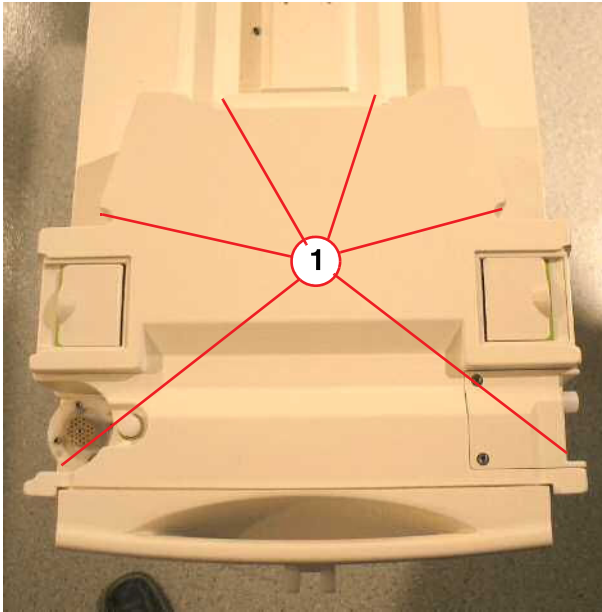

(1) Screws to be removed


- Remove the carrier plate of the RTT.
- Remove six screws on top of the cover from the coil plugs
- Disconnect cables.

Fig. 116: Foot-end cover



- Remove the two M8 Allen screws, if necessary disconnect cables, and remove foot-end plate.

Fig. 117: Mechanical end stop

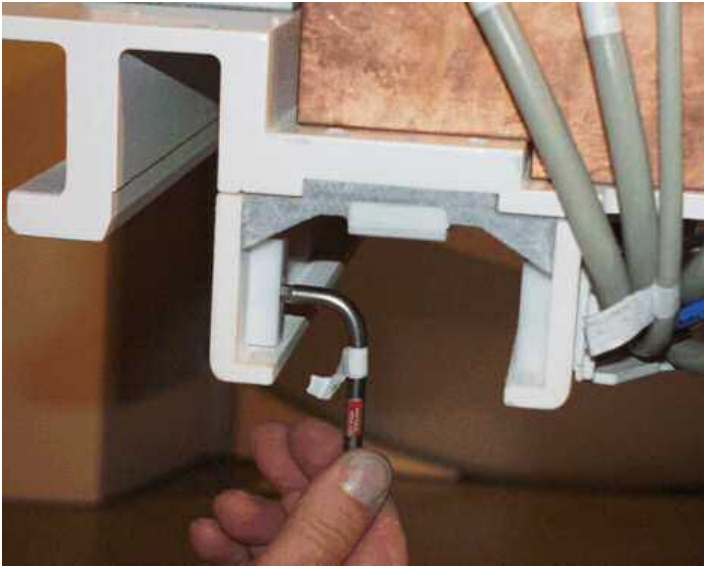

- Remove mechanical end stops.

3.5.4.2 Cables/energy chain

- Disconnect all cables and air hoses. Use a flat screwdriver to open the RF connectors.

Fig. 118: Cable connections



Fig. 119: Energy chain



- Dismount the energy chain.

3.5.4.3 Tabletop plate

- Release the emergency clutch and push the tabletop plate manually into the magnet until the cantilever arm is no longer touching the tabletop plate.

Fig. 120: Tabletop plate, pushed into the magnet


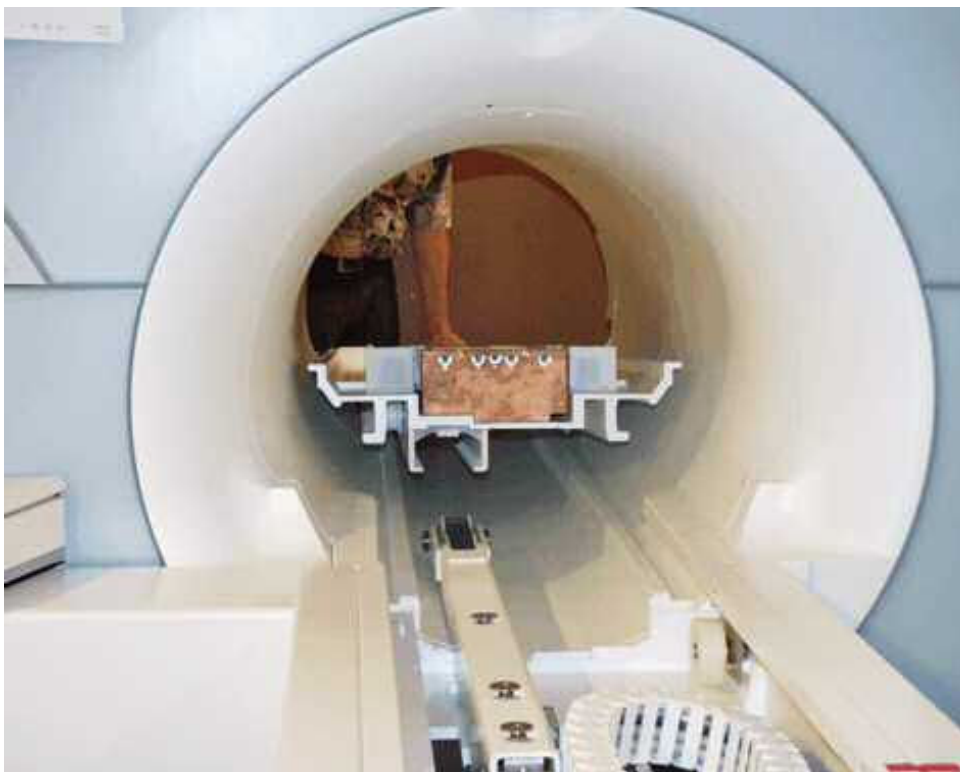


Fig. 121: Tabletop moved into the magnet



- Carefully push the tabletop plate manually into the magnet; not valid for Trio A Tim Systems.

Fig. 122: Tabletop, removed

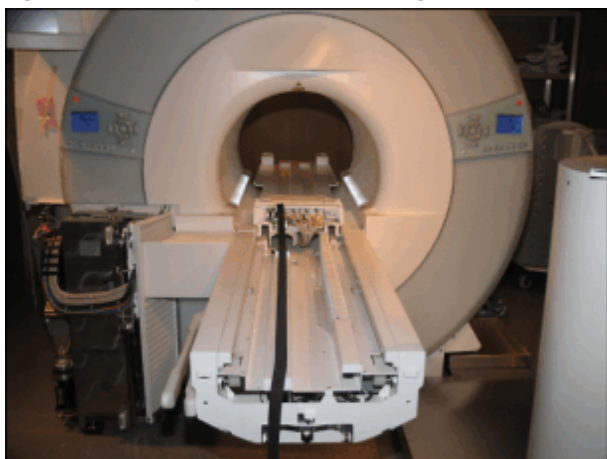

- Lift up the tabletop - a second person is necessary - and carry it to the front of the magnet, replace tabletop plate.


3.5.5 Assembly

For assembling the partial drive, please refer to (→ Assembly / Page 167), not valid for Trio A Tim Systems.




Be careful to avoid pinching any cables when assembling the cover of the coil plug connectors.

Fig. 123: Distance of cantilever arm

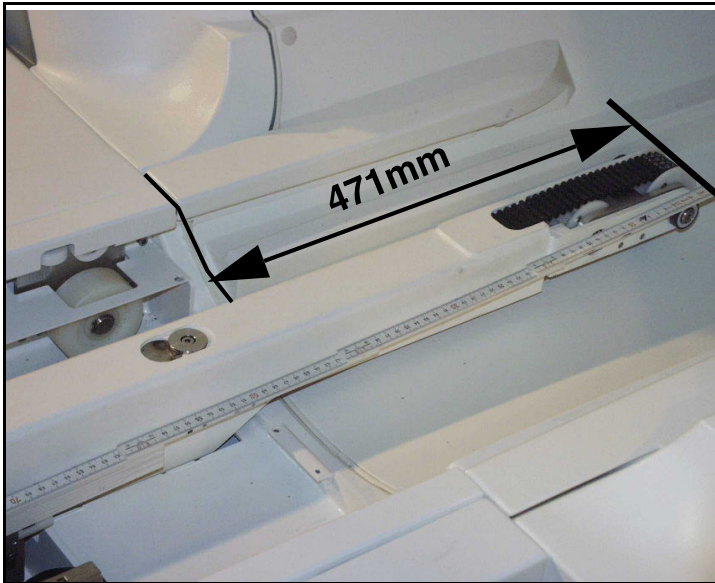

- Before you insert the new tabletop plate: Please do a check of the distance from the end of the cantilever arm to the support frame and adjust to 471 mm.
- Move the new tabletop plate from the front to the back of the magnet.
- Install the energy chain.

Fig. 124: Mark on the telescopic arm

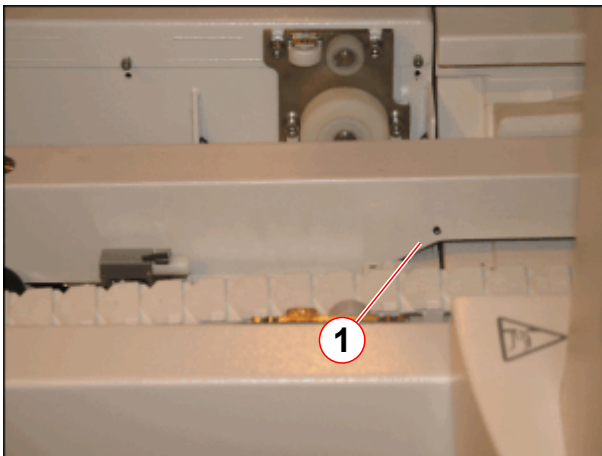

(1) Marking hole in the telescopic arm

- On newer systems the distance mark shows the correct length of the cantilever arm.



Use the same length of screws when remounting the energy chain.

3.5.6 Concluding steps

- Install the end stops according to the system/room specification, (Trio, A Tim/Symphony, A Tim System, may reduce the movement range of the horizontal movement).
- Reconnect all cables and air hoses.



Be careful to avoid pinching any cables when assembling the covers of the patient table.

- Install foot-end plate and the cover of the coil plugs.

3.5.7 Adjustments

For adjustment of the toothed belt, please refer to (→ Adjustments / Page 170)

3.5.8 Final checks

- Switch ON F22 (RF-room) ptab and check for smooth horizontal movement of the tabletop plate.

3.6 Body of the energy chain

3.6.1 General

All cables of the ptab are running in a flexible chain (energy chain). This chain can be replaced.

3.6.2 Tools and time required

3.6.2.1 Time

Approx. 30/45 min

3.6.2.2 Tools

Non-magnetic tool kit

3.6.3 Preliminary

- Switch OFF F22 (RF-room) in the ACC.



According to the system/room specifications the movement range of the horizontal movement can be reduced

With the whole body version of the horizontal drive, the movement range of the horizontal movement may be reduced. If that is the case, remove the end stops.
(→ Fig. 289 Page 173)

Before disassembly the patient table covers on a system with a partial drive unit: The toothed belt must be released to push the table top to the service end of the magnet. Therefore refer to. (→ Prerequisite for cover removal / partial drive / Page 98)

3.6.4 Disassembly

3.6.4.1 Covers



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

Basic table

Fig. 125: Coil plug cover (foot end)



- Remove the Allen screws on top of the cover from the coil plugs.


Fig. 126: Remove Allen screws



- Remove the lower Allen screws of the foot-end cover.

Removable table


Fig. 127: RF-Cover RTT

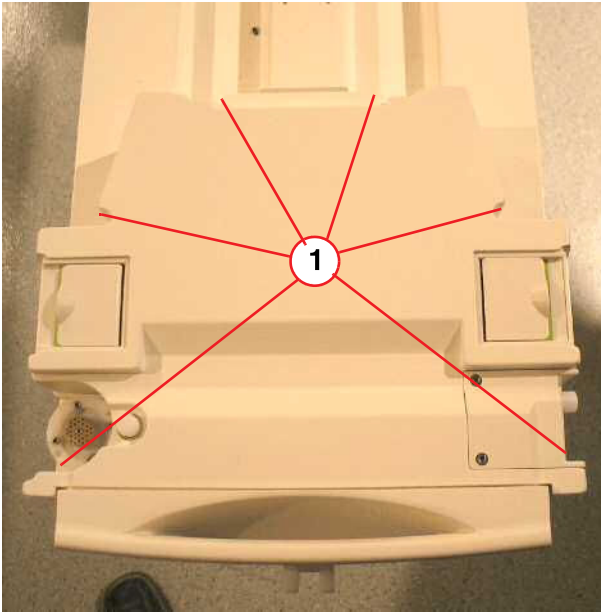

(1) Screws to be removed

- Remove the carrier plate of the RTT.
- Remove six screws on top of the cover from the coil plugs.
- Disconnect cables.

Fig. 128: Remove end plate

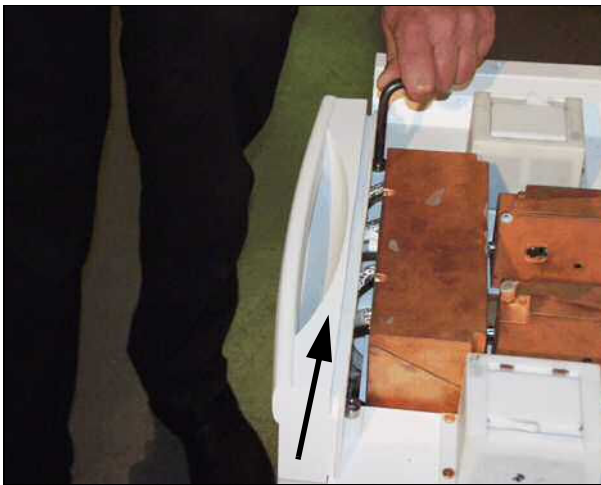

- Remove the two M8 Allen screws, if necessary disconnect cables, and remove foot-end plate.


3.6.4.2 Energy chain

Fig. 129: Energy chain



- Push the tabletop plate manually into the magnet to access the energy chain.
- Remove the lower part of the energy chain.

Fig. 130: Remove cover of energy chain

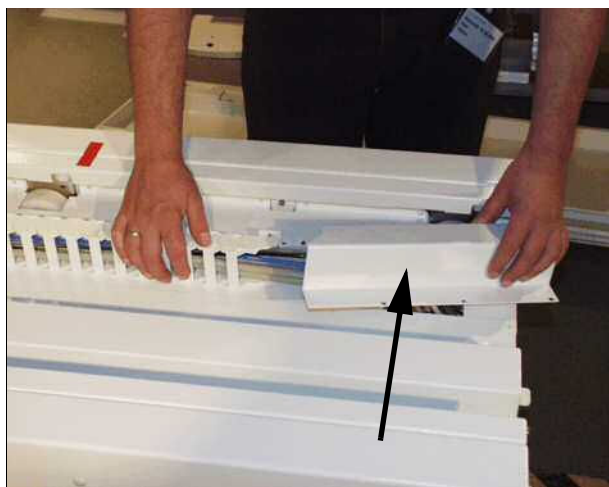

- Remove the upper energy chain cover.

Fig. 131: Remove energy chain



- Remove the upper energy chain of the bracket.

3.6.5 Assembly




Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

- Open each chain link by using a screwdriver and replace damaged parts



Make sure all cables are routed in parallel inside the chain

Fig. 132: Release the chain links

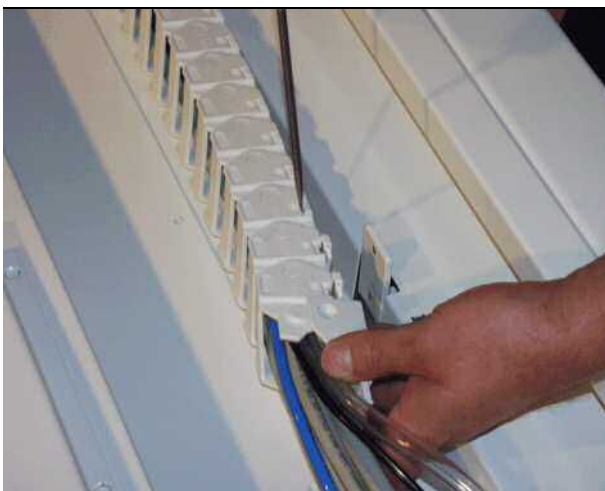

- Install the energy chain.

3.6.6 Concluding steps

- Install covers.

3.6.7 Adjustments

For adjustment of the toothed belt, please refer to (→ Adjustments / Page 170), not valid for Trio A Tim Systems.

3.6.8 Final checks

- Switch ON F22 (RF-room) and check for smooth movement of the energy chain

3.7 Guide rollers

3.7.1 General

The table is moved in the horizontal direction via guide rollers (special bearings). Each roller can be replaced separately.

3.7.2 Tools and time required

3.7.2.1 Time

Approx. 30/45 min

3.7.2.2 Tools

Non-magnetic tool kit
Loctite 221

3.7.3 Preliminary

- Switch OFF F22 (RF-room) in the ACC



According to the system/room specifications the movement range of the horizontal movement can be reduced.

With the whole body version of the horizontal drive, the movement range of the horizontal movement may be reduced. If that is the case, remove the end stops.
(→ Fig. 289 Page 173)

Before disassembly the patient table covers on a system with a partial drive unit: The toothed belt must be released to push the table top to the service end of the magnet. Refer to (→ Prerequisite for cover removal / partial drive / Page 98)

3.7.4 Disassembly



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

3.7.4.1 Covers

- Push the tabletop plate manually into the magnet to access the guide rollers.

- Remove Allen screws and remove the affected side covers.

Fig. 133: Remove cover from the affected side

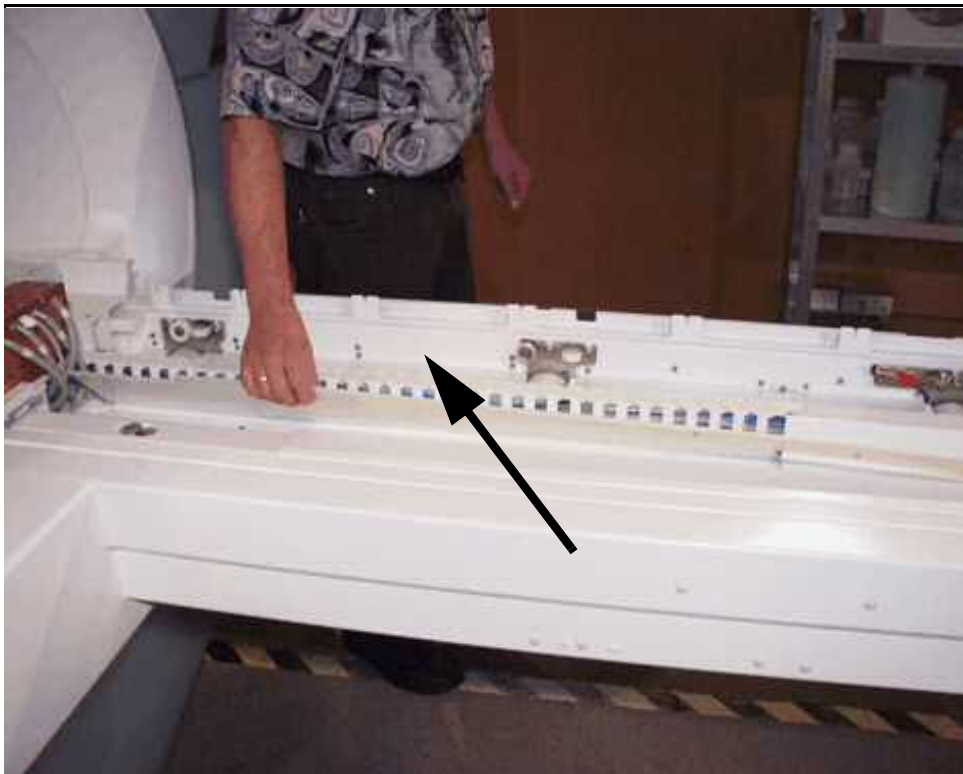


Fig. 134: Remove cover from the affected side


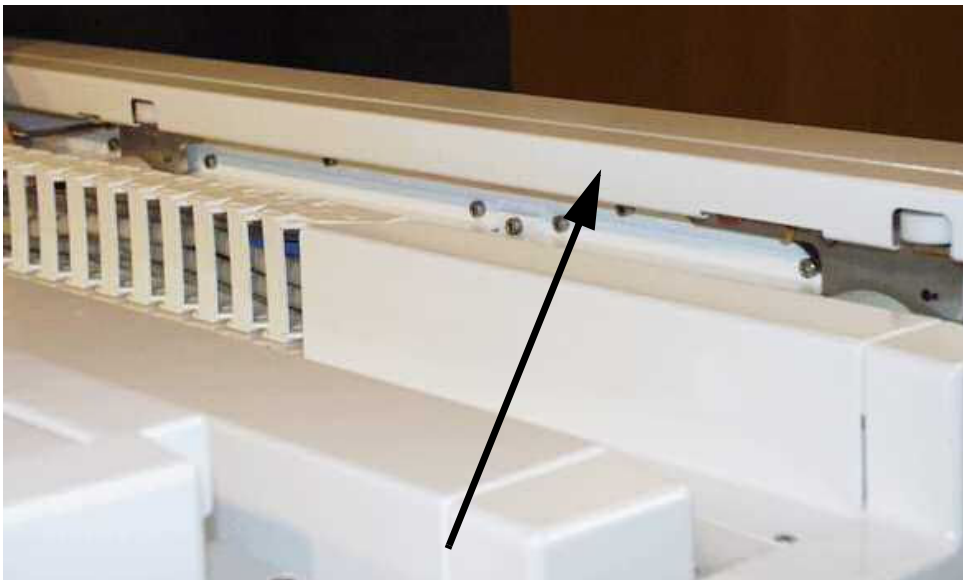


Fig. 135: Remove cover of energy chain




- Remove the upper energy chain cover.

Fig. 136: Remove energy chain

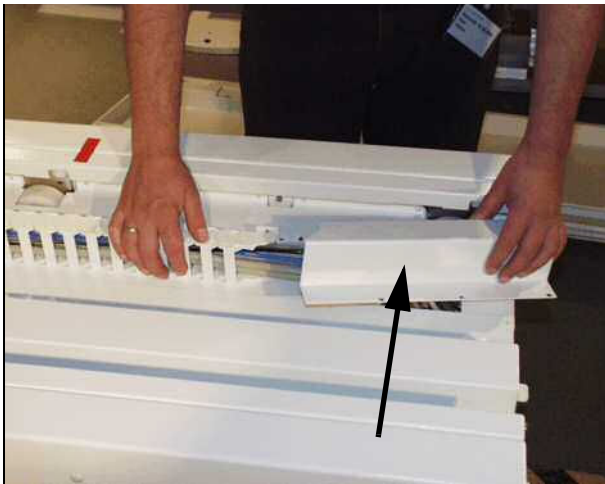

- Remove the upper energy chain of the bracket.

Fig. 137: Foot-end cover

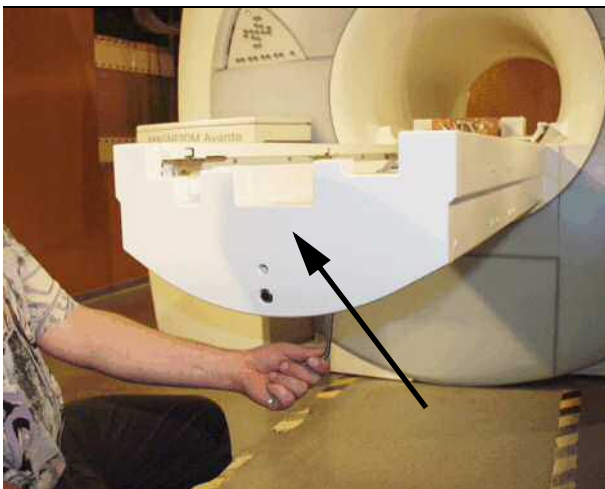


- Remove the foot-end cover of the patient table.

3.7.5 Assembly

For assembling the partial drive, please refer to (→ Assembly / Page 167) not valid for Trio A Tim Systems.

- Replace the defective guide roller and safety screws by apply Loctite 221.

Fig. 138: Replace affected guide roller

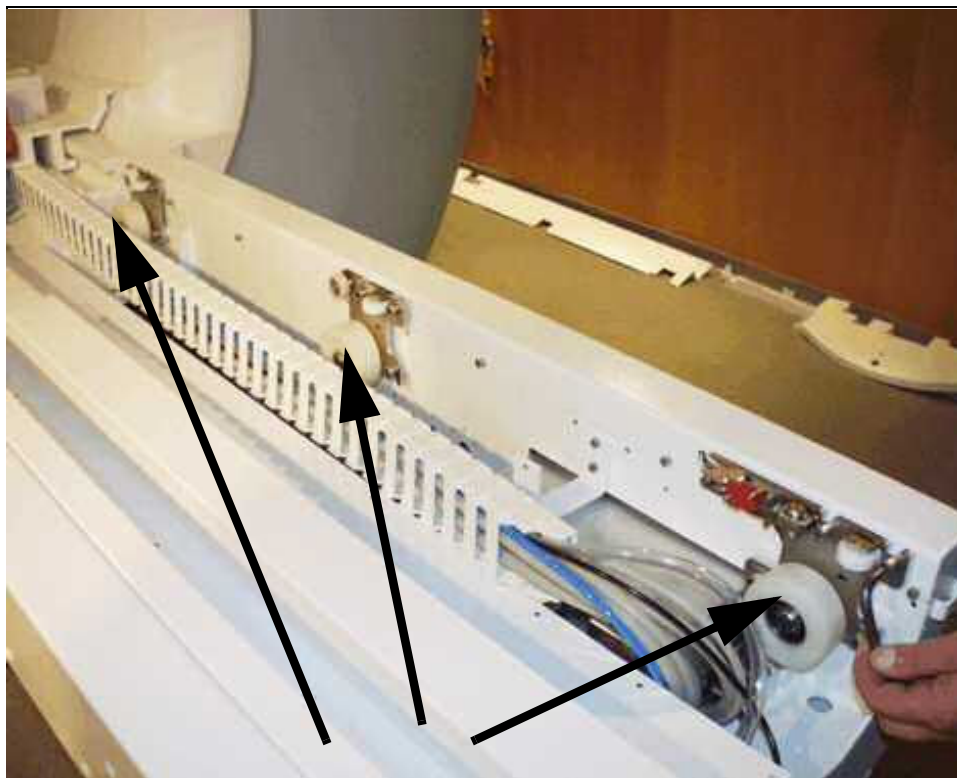

3.7.6 Concluding steps

- Switch ON F22 (RF-room) in the ACC.

3.7.7 Adjustments

For adjustment of the toothed belt, please refer to (→ Adjustments / Page 170), not valid for Trio A Tim Systems.

3.7.8 Final checks

- Install covers.
- Check smooth horizontal movement.

3.8 MDSD device

3.8.1 General



CAUTION

Strong magnetic field!

Magnetic objects may cause injury.

- » **Hold the ferromagnetic part with both hands while moving it to the correct position and to remove the unit directly out of the field!**

The MDSD device is responsible for driving the horizontal movement of the ptab. Two types of MDSD are introduced.

MDSD for whole body movement.



With revision level 2 and higher, the grounding shield and the external filter must not be used anymore. On older systems, the grounding shield and/or the external filter may be installed. Partial drive is not affected.

MDSD for partial drive, not valid for Trio A Tim Systems.

Each of the MDSD device can be replaced, **they are n o t compatible.**

3.8.2 Tools and time required

3.8.2.1 Time

Approx. 120 min


3.8.2.2 Tools

Non-magnetic tool kit

Philips screwdriver

Service aid MDSD (delivered with the MDSD spare part)

Fig. 139: MDSD service aid

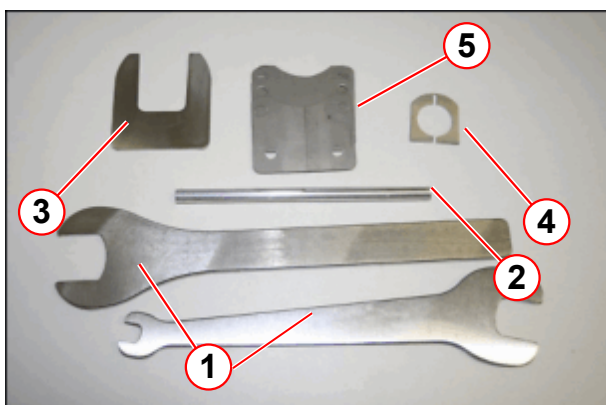

- (1) Spanners
- (2) Feed-in tool for fiber optic cable
- (3) Clamp claw
- (4) Clamp collar
- (5) Shim plates



For patient tables with partial drive MDSD A4330 (p.n. 10277315), the only service tool included in the spare part is the feed-in tool fiber optical cable. Only for the whole body MDSD A4320 (p.n. 8112240) a complete service aid set, as shown in the figure above, is needed and included in the spare part.

3.8.3 Preliminary

- Position patient table at home position
- Switch OFF F22 (RF-room) in the ACC



According to the system/room specifications the movement range of the horizontal movement can be reduced.

With the whole body version of the horizontal drive, the movement range of the horizontal movement may be reduced. If that is the case, remove the end stops.
(→ Fig. 289 Page 173)

Before disassembly the patient table covers on a system with a partial drive unit: The toothed belt must be released to push the table top to the service end of the magnet. Therefore refer to (→ Prerequisite for cover removal / partial drive / Page 98), not valid for Trio A Tim Systems

3.8.4 Disassembly

3.8.4.1 Prerequisite for cover removal / partial drive




Be careful and avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.



The tabletop plate may fall - support the tabletop plate.

Basic table

Fig. 140: Coil plug cover (foot end)

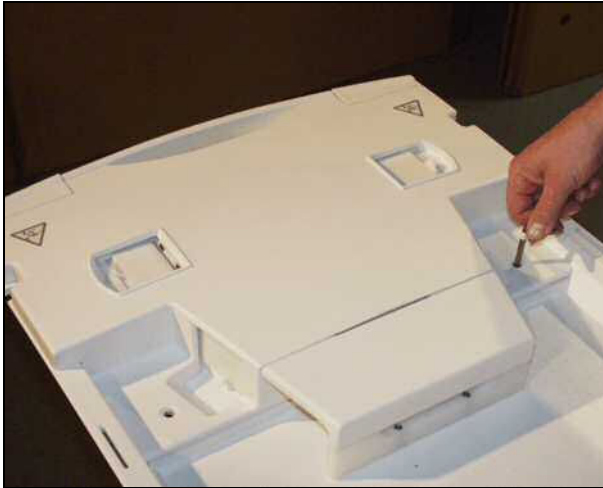

- Remove the Allen screws on top of the cover from the coil plugs.


Fig. 141: Remove Allen screws



- Remove the lower Allen screws of the cover from the coil plugs.

Removable table


Fig. 142: RF-Cover RTT

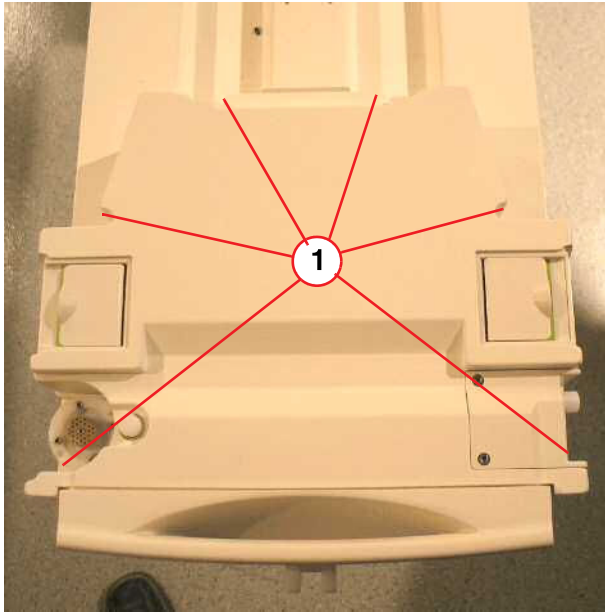

(1) Screws to be removed

- Remove the carrier plate of the RTT.
- Remove six screws on top of the cover from the coil plugs.
- Disconnect cables.

Fig. 143: Additional mounting bracket for foot end cover (newer systems)

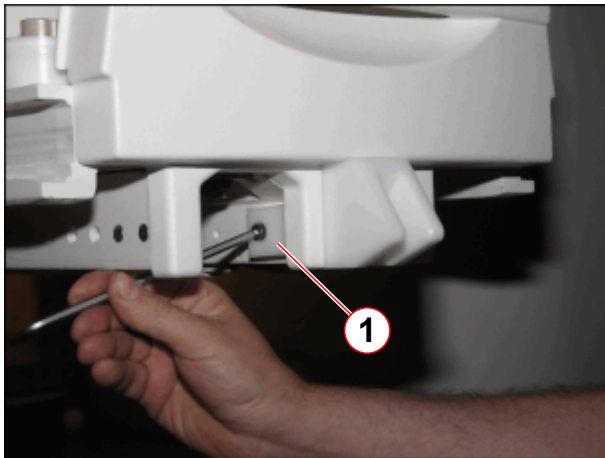

(1) Allen screw


- On newer systems, remove Allen screw of mounting bracket.

Fig. 144: Foot-end cover



- Remove screws of the foot-end cover.


Fig. 145: Foot end cover, removed

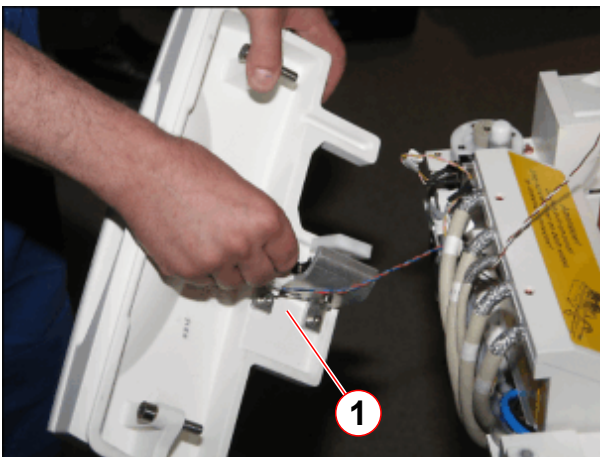

- Disconnect cable connections and remove the foot-end cover..
- Slide out the lower cover on the bottom of the patient table.

(1) Disconnect cables


Fig. 146: Mechanical end stop of partial drive unit

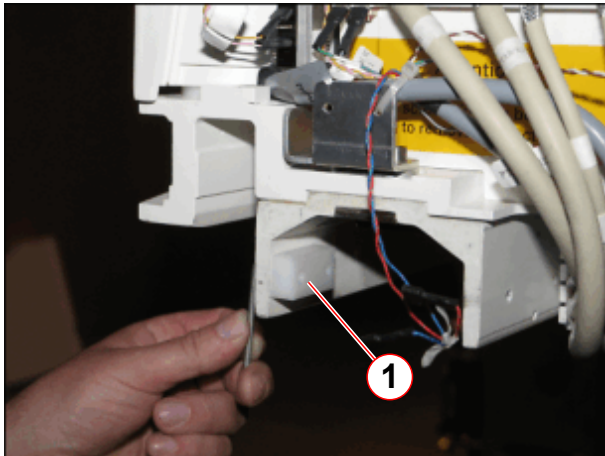

(1) End stop to be removed

- Remove mechanical end stop.


Fig. 147: Tension unit for toothed belt

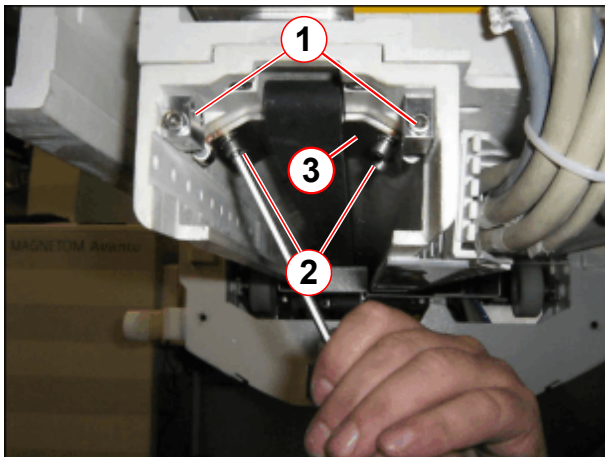

(1) Allen screws for adjusting belt tension
(2) Attachment screws
(3) clamping plate

- Release tension of toothed belt (1), remove Allen screws (2) and tooth belt tension unit, not valid for Trio A Tim Systems.


Fig. 148: Cover of toothed belt mount (head end)

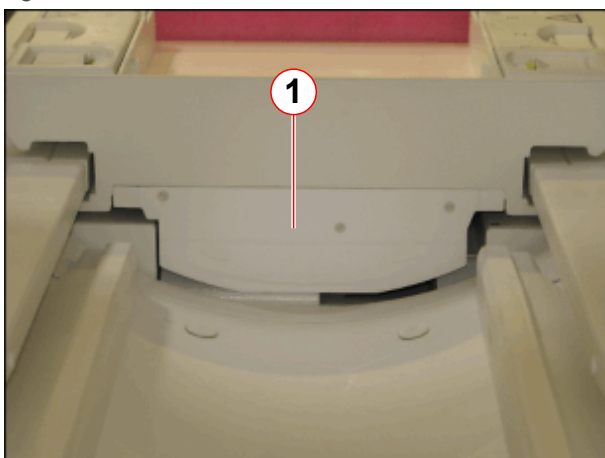

(1) Cover of toothed belt mount

- Remove the cover of the toothed belt mount of toothed belt (1) on the head end from the ptab, not valid for Trio A Tim Systems.


Fig. 149: Toothed belt attachment (head end)

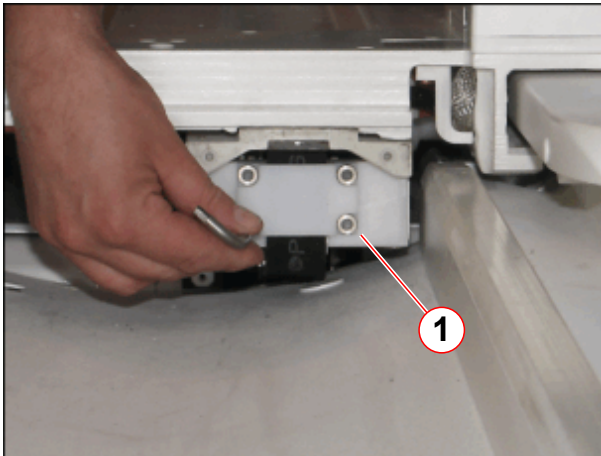

(1) Allen screws of toothed belt mount

- Remove the toothed belt mount (1) on the head end from the ptab, not valid for Trio A Tim Systems.

Fig. 150: Tabletop moved into the magnet

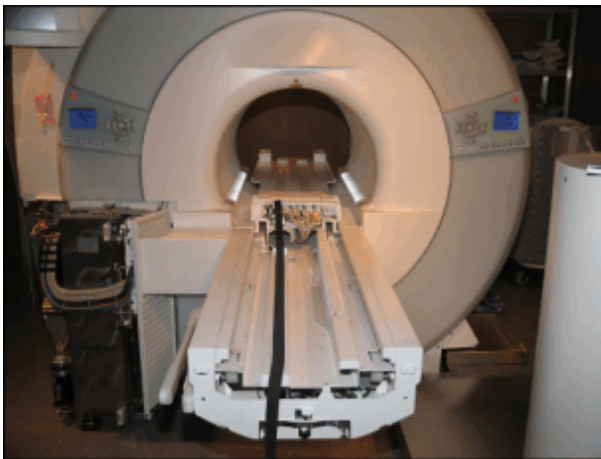

- Carefully push the tabletop plate manually into the magnet.

3.8.4.2 Covers




Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.



For disassembling/assembling the patient table covers we recommend a mid-sized Philips screwdriver.


- Push the tabletop plate manually into the magnet (whole body movement).

Fig. 151: Foot-end cover

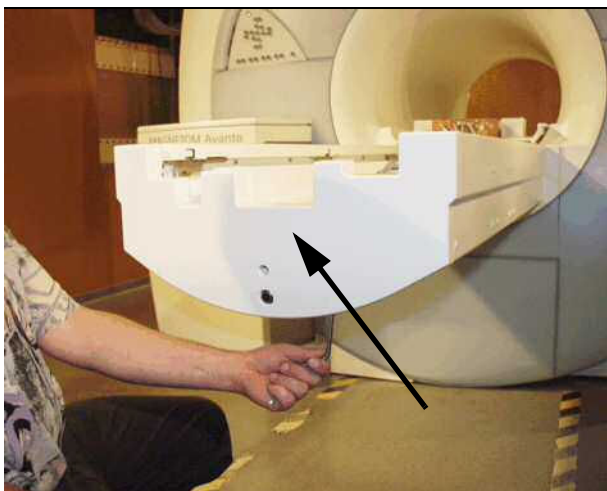

- Remove the foot-end cover of the patient table.


Fig. 152: Plastic inner cover

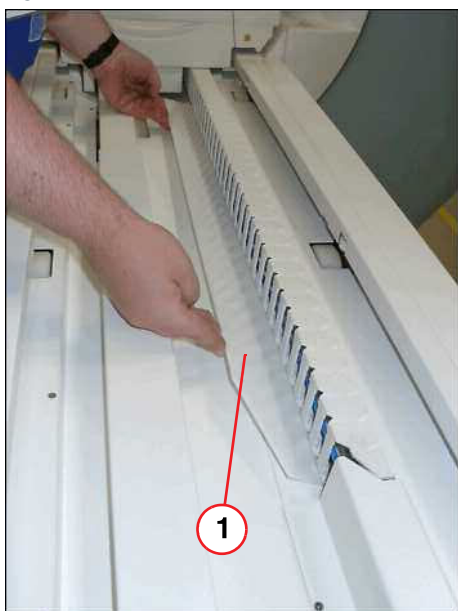

- Remove inner plastic cover.

(1) Plastic inner cover

Fig. 153: Remove cover of energy chain

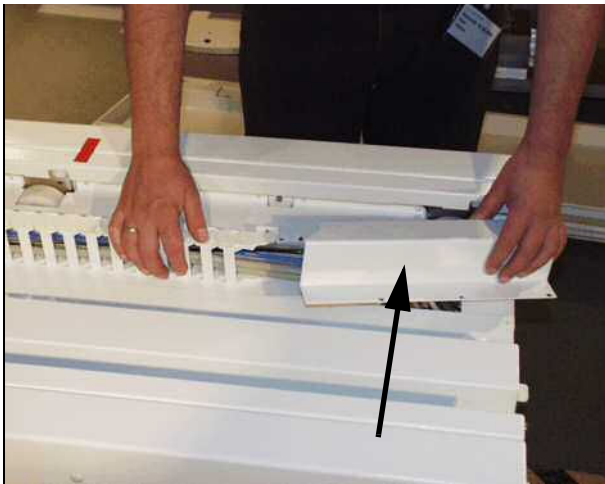

- Remove the upper energy chain cover.


Fig. 154: Remove energy chain



- Remove the upper energy chain of the bracket.

- Remove right guide rail cover (do not remove screws, only loosen).


Fig. 155: Inner main cover

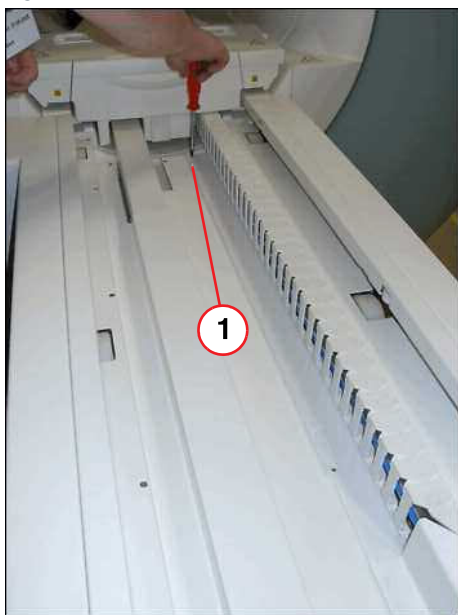

(1) Screws to be removed

- Remove screws of main inner cover.


Fig. 156: Inner cover

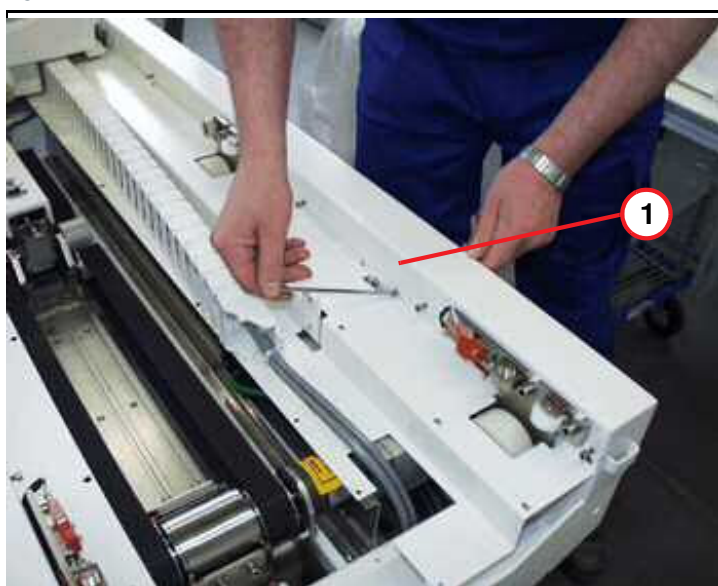

(1) Top inner cover

- Remove right inner cover.

Fig. 157: Allen screws

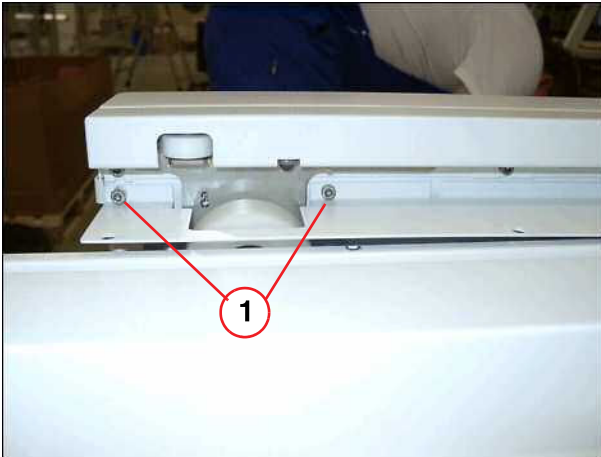


(1) Screws to be removed

- Remove two Allen screws on the left side on the table end (foot) before you remove the main inner cover.


- Remove bracket above the MDSD device.
- Remove main inner cover.

Fig. 158: Disassembling

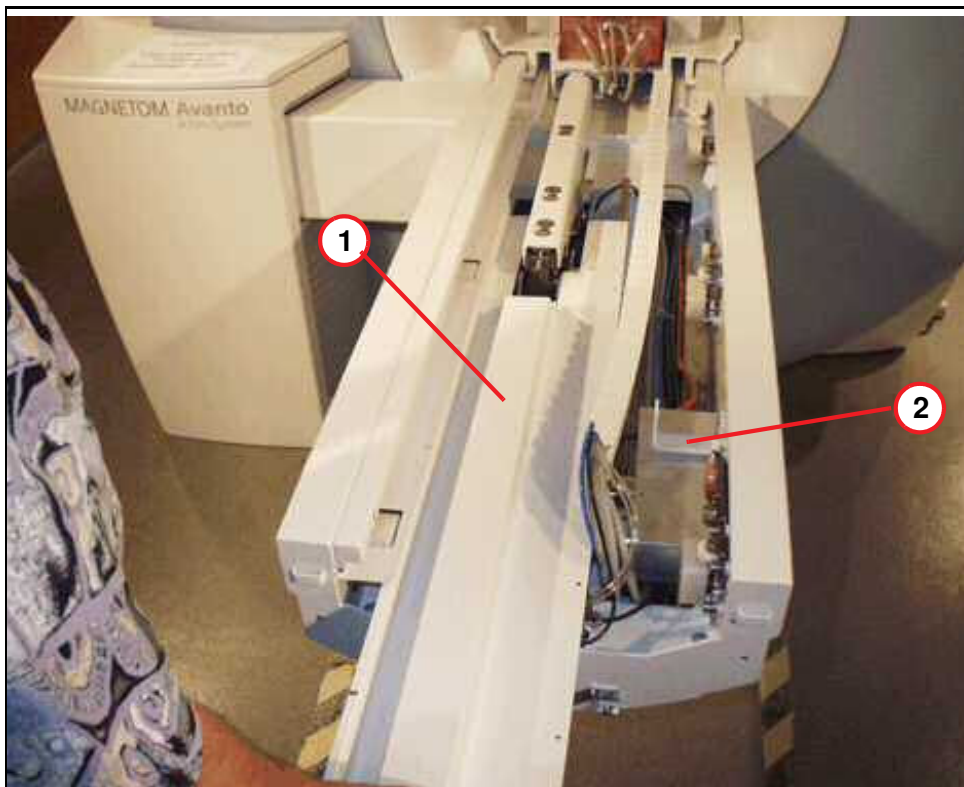

- (1) Main inner cover
(2) Mounting bracket for energy chain

Fig. 159: Inner cover, left side

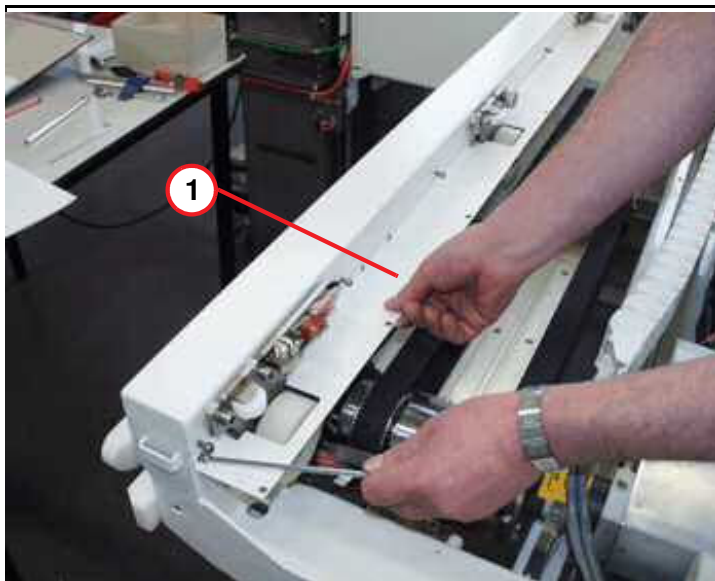


(1) Inner cover, left side

- Covers removed.

3.8.4.3 MDSD device of whole body drive


Fig. 160: Example - removing FOCs from MDSD

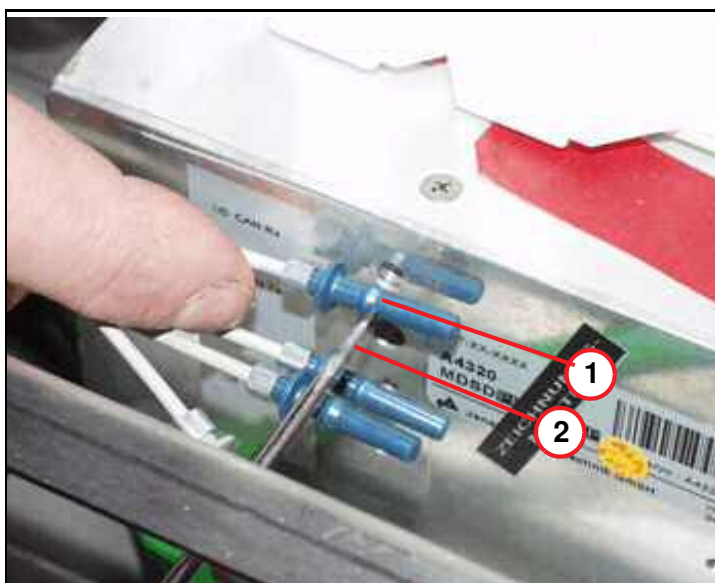

(1) Cam of FOC connector
(2) Screwdriver, flat blade

- The figure below shows the FOC plug connection (principle) inside the MDSD device. Use a flat blade screwdriver to disconnect the FOCs inside the MDSD device.

Fig. 161: Fiber-optic cables of MDSD

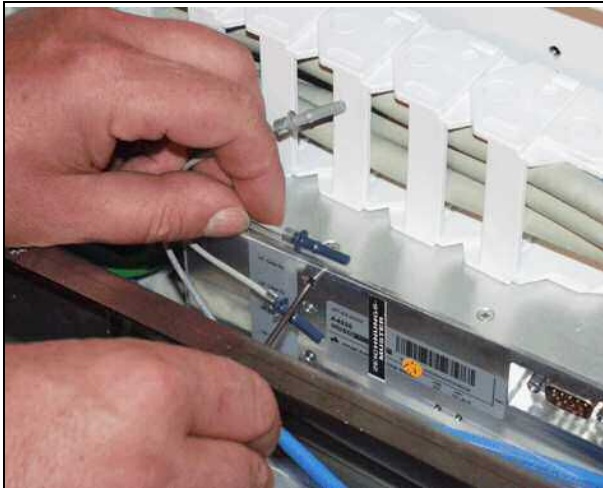


- Feed the flat blade screwdriver through the hole in the aluminum cover and into the FOC's plug connection. Try to reach the cam of the FOC and carefully push it sideways while at the same time pulling it out from the FOC. Refer to illustration.

- Release the horizontal emergency movement, enabling manual horizontal movement.




To avoid the overwinding of the MDSD clamping nut, the clamp collar and the claw fastener (part of the service aid) must be used as described.

Fig. 162: MDSD service aid, principle of use

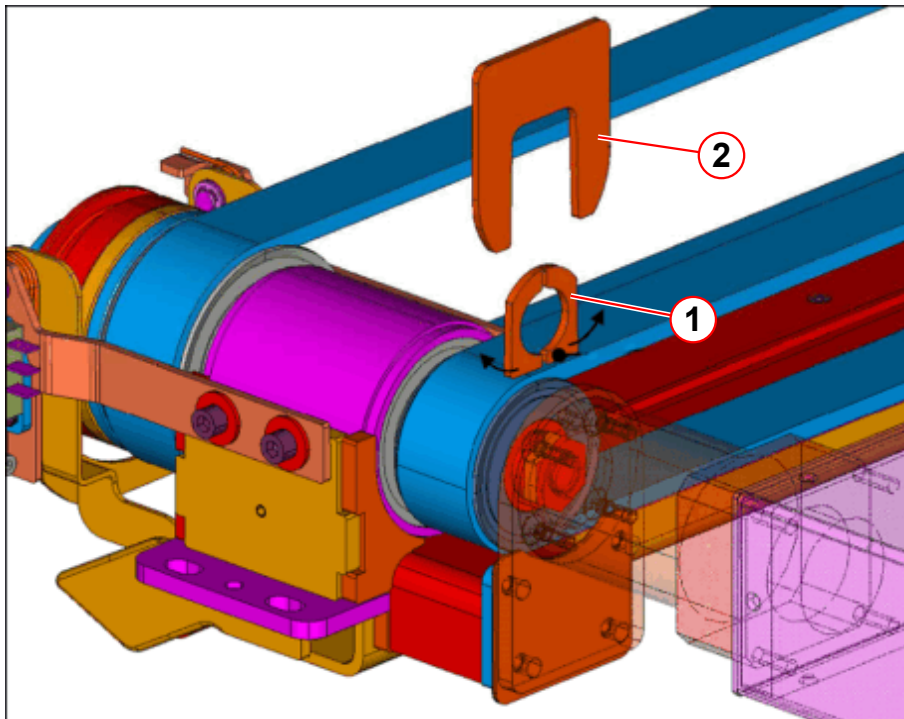


- (1) Clamp collar
(2) Claw fastener

Fig. 163: Clamp collar

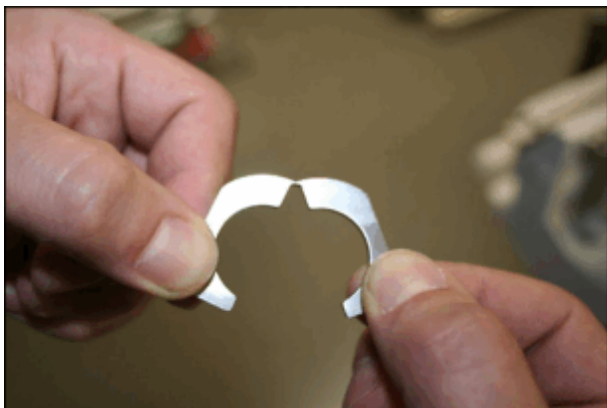

- Carefully bend up the clamp collar, as shown.

Fig. 164: Clamp collar

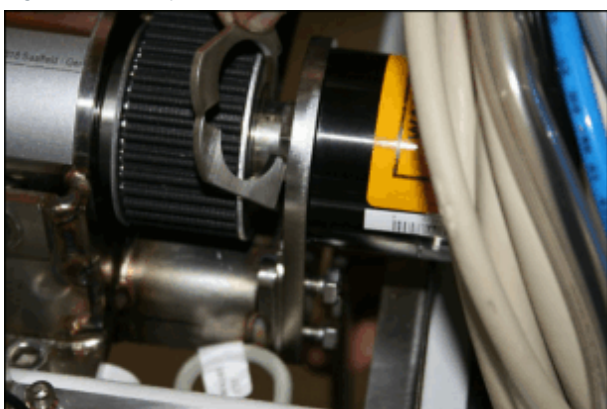

- Insert the clamp collar.

Fig. 165: Clamp collar



- Press the clamp collar together.


(1) Insert clamp collar

Fig. 166: Claw fastener



(1) Insert claw fastener over the clamp collar.

Fig. 167: RELEASE or TIGHTENING of the clamping nut

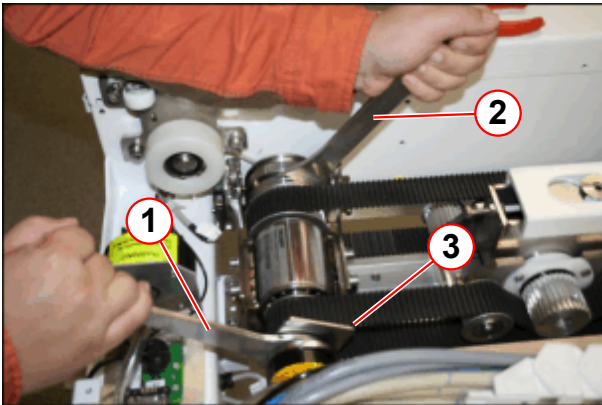

- (1) RELEASE, pulling the spanner toward the foot end, TIGHTEN, pushing toward the magnet
- (2) Holding the spanner tight
- (3) Claw fastener

- Insert the claw fastener (1) over the clamp collar.
- Use the open-end spanners of the MDSD service aid to release the MDSD clamping nut.
- **RELEASE**, pulling the open-end spanner (1) **toward the foot end**, while holding the open-end spanner (2).
- The claw fastener should remain.

Fig. 168: *RELEASING or TIGHTENING the clamping nut*


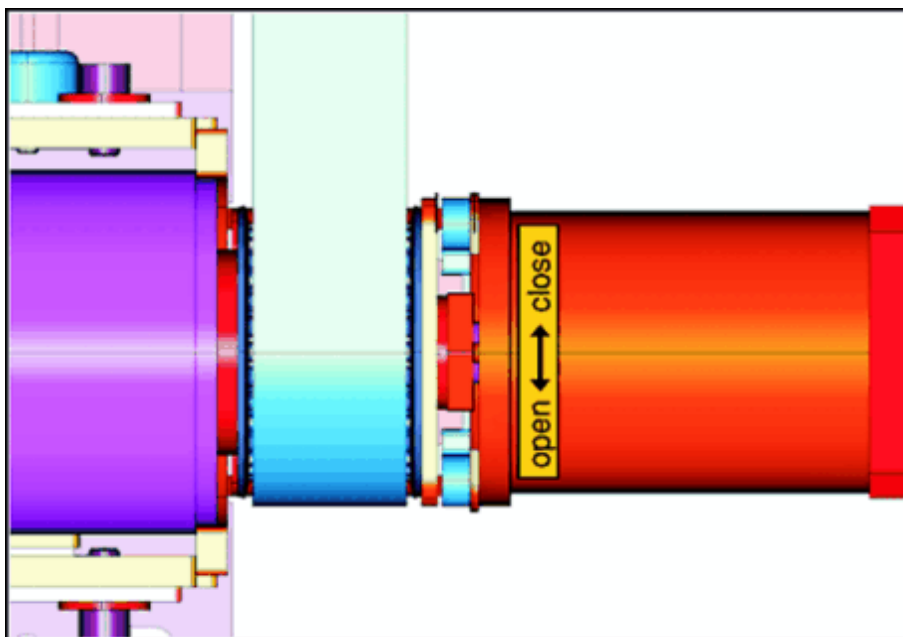



Fig. 169: *M6 hex screws of MDSD bracket*

- Remove the four M6 hex nuts.

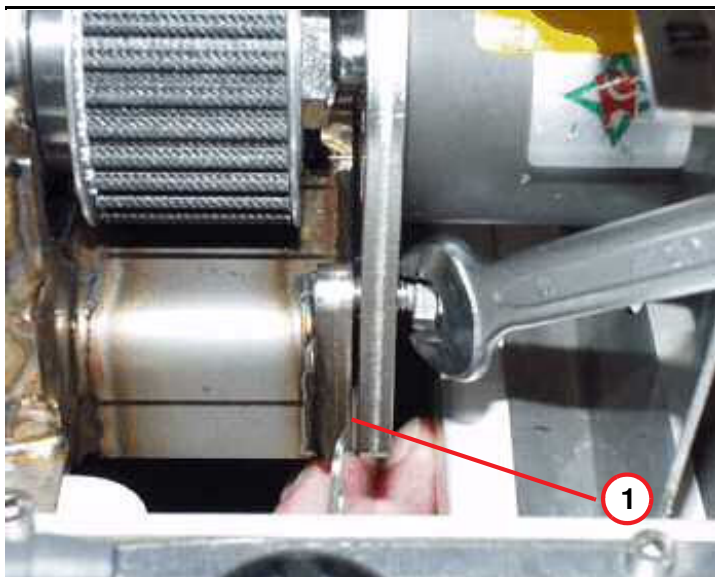

(1) Shim gap appears when securing the 22 mm nut

Fig. 170: MDSD device

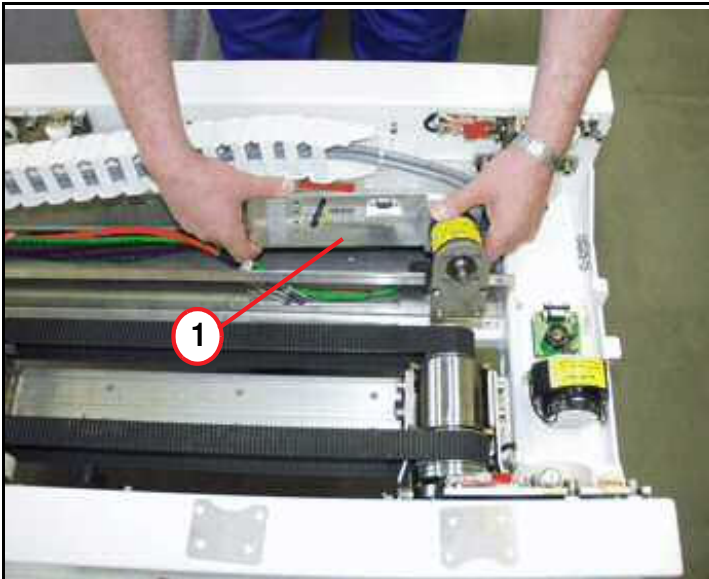

(1) Replace MDSD device

- Replace the MDSD drive, take care of the correct version.
- Remove the MDSD service aid (clamp collar and claw fastener).

! CAUTION

The MDSD drive consists of magnetic objects!
Magnetic objects may cause injury.


- » Maintain a safe distance to the magnetic field!

! CAUTION

Strong magnetic field!
Magnetic objects may cause injury.

- » Hold the ferromagnetic part with both hands while moving it to the correct position and to remove the unit directly out of the field!

Fig. 171: MDSD with external filter and grounding shield

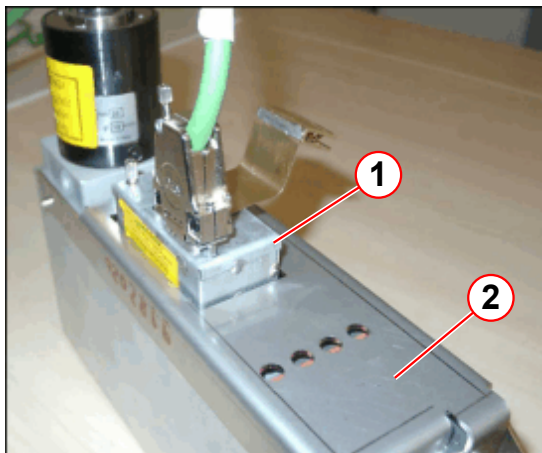


- (1) External filter
- (2) Grounding shield

Note:

If a grounding shield and/or an external filter is installed, **do not use them anymore for the new MDSD** (applies to revision level 02 and higher).

3.8.4.4 MDSD device of partial drive, not valid for Trio A Tim Systems


Fig. 172: Clutch of MDSD, partial drive

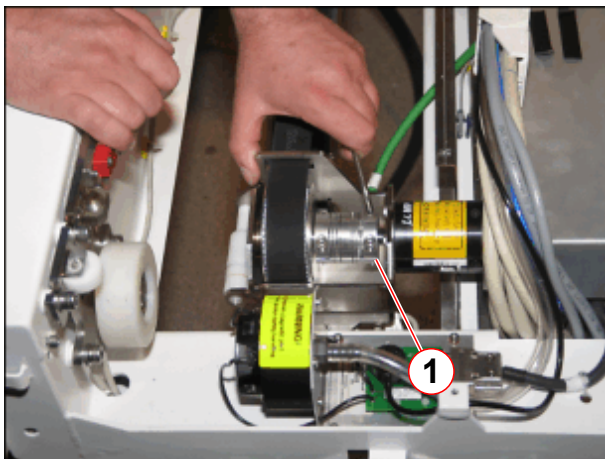

- (1) Allen screws of MDSD clutch


- Remove Allen screws of MDSD clutch (1).

Fig. 173: Clutch of MDSD removed

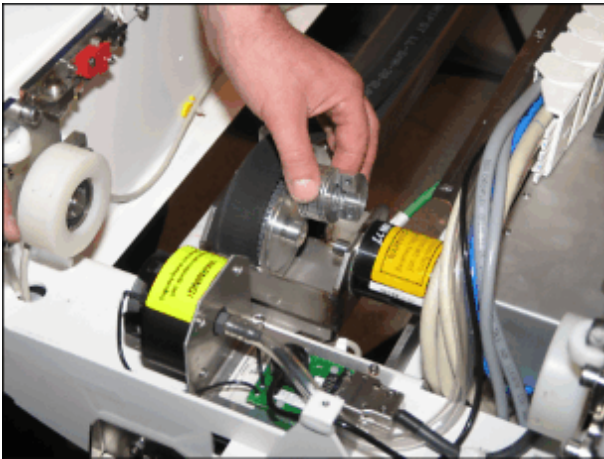

- Remove MDSD clutch.

Fig. 174: MDSD clutch removed

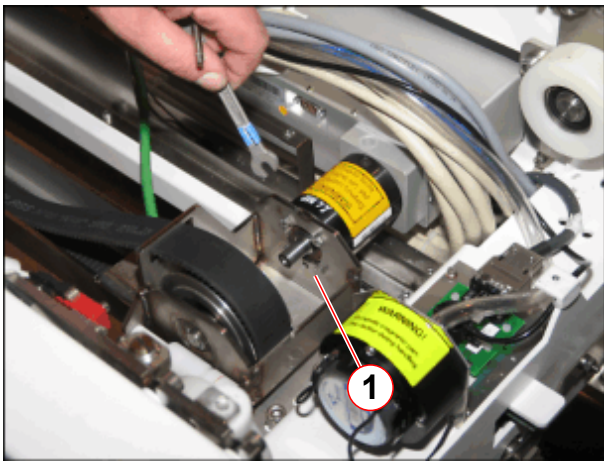

(1) Mounting screws of MDSD

- Remove the hex screws (1) of the mounting bracket from the MDSD.

Fig. 175: Removing the MDSD



- Replace the MDSD drive.
- Hold the ferromagnetic part with both hands while moving it to the correct position and to remove the unit directly out of the field!

**CAUTION**

The MDSD drive consists of magnetic objects!

Magnetic objects may cause injury.

- » Maintain a safe distance to the magnetic field.

**CAUTION**

Strong magnetic field!

Magnetic objects may cause injury.

- » Hold the ferromagnetic part with both hands while moving it to the correct position and to remove the unit directly out of the field!

3.8.5 Assembly

3.8.5.1 MDSD drive of whole body drive




The service aid for the MDSD is a prerequisite for assembly.



To avoid the overwinding of the MDSD clamping nut, the clamp collar and the claw fastener (part of the service aid) must be used as described.

Fig. 176: MDSD service aid, principle of use

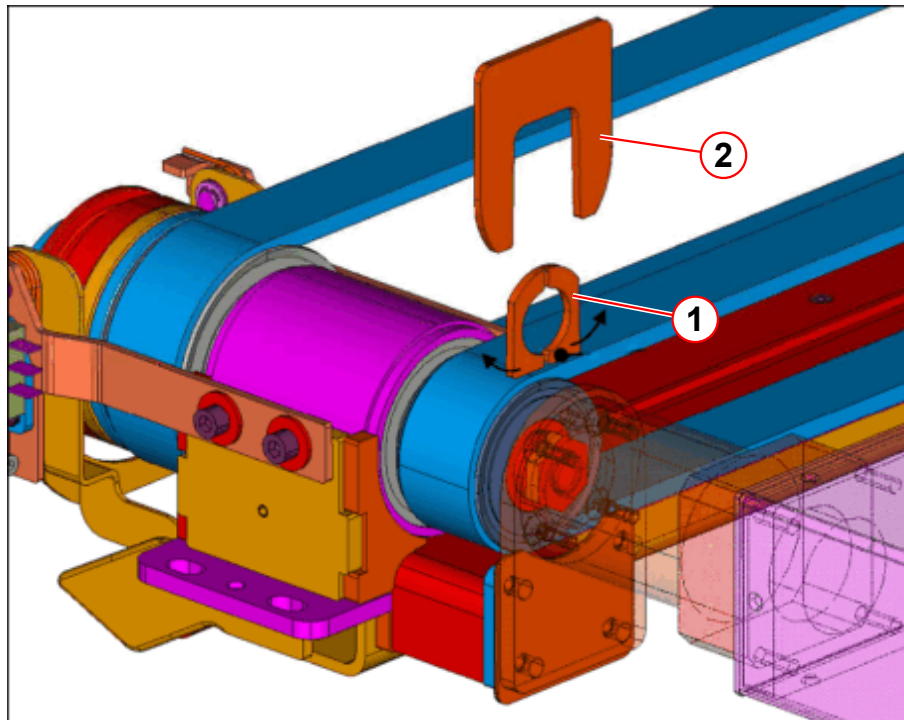


- (1) Clamp collar
- (2) Claw fastener

Fig. 177: Set of shim plates MDSD drive



- Shim plates must be used when assembling the MDSD device.


Fig. 178: Hex screws of MDSD

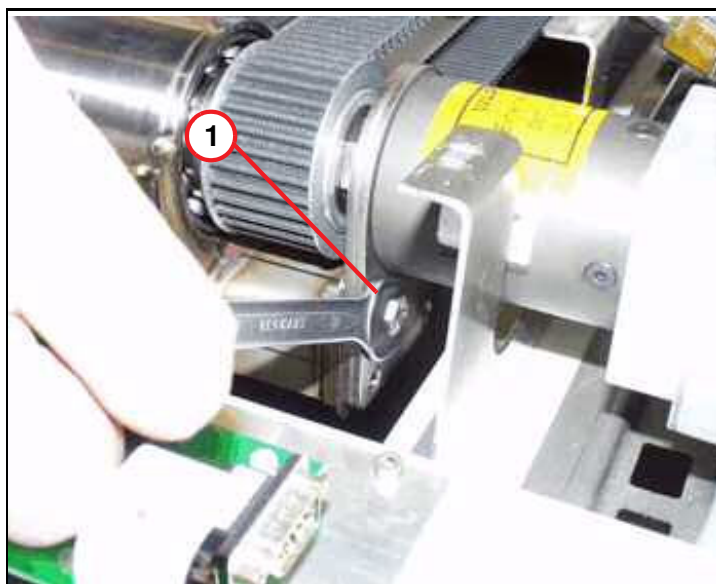

(1) M6 hex screws


- Insert only the upper two M6 hex nuts and don't tighten them, use them only to hold (as a torque bracket).

Fig. 179: Clamp collar

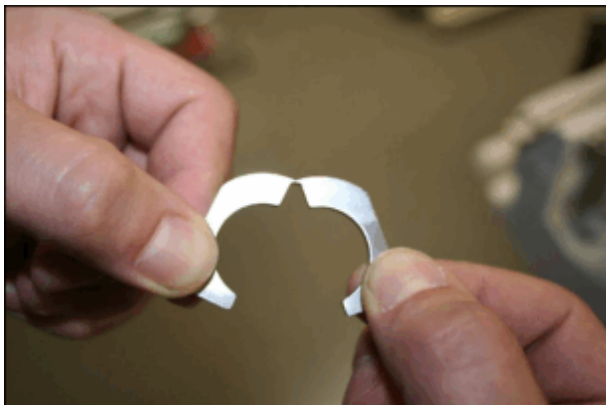


- Carefully bend up the clamp collar, as shown.

Fig. 180: Clamp collar

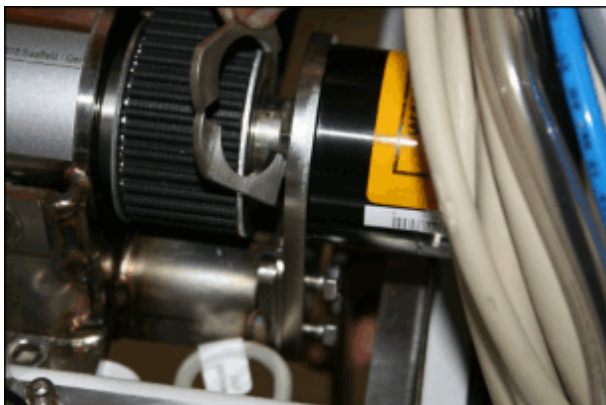

- Insert the clamp collar.

Fig. 181: Clamp collar

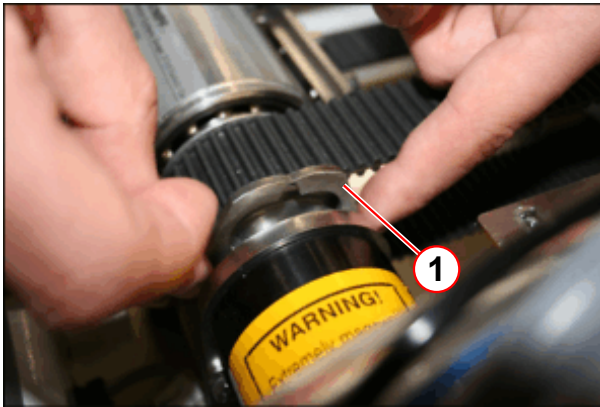


(1) Insert clamp collar

Fig. 182: Claw fastener



(1) Insert claw fastener over the clamp collar.

Fig. 183: RELEASE or TIGHTENING of the clamping nut

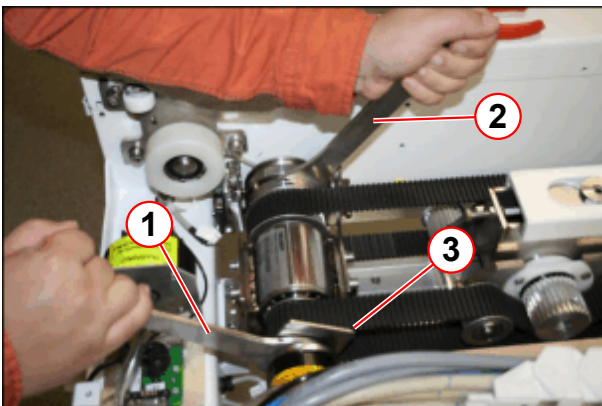

- (1) RELEASE, pulling the spanner toward the foot end, TIGHTEN, pushing toward the magnet
- (2) Holding the spanner tight
- (3) Claw fastener

- Press the clamp collar together.
- Insert the claw fastener (1) over the clamp collar.
- Use the open-end spanners of the MDSD service aid to tighten the MDSD clamping nut.
- **TIGHTEN**, pushing the open-end spanner (1) **to the magnet**, while holding the open-end spanner (2).
- The claw fastener should remain.

Fig. 184: *RELEASING or TIGHTENING the clamping nut*


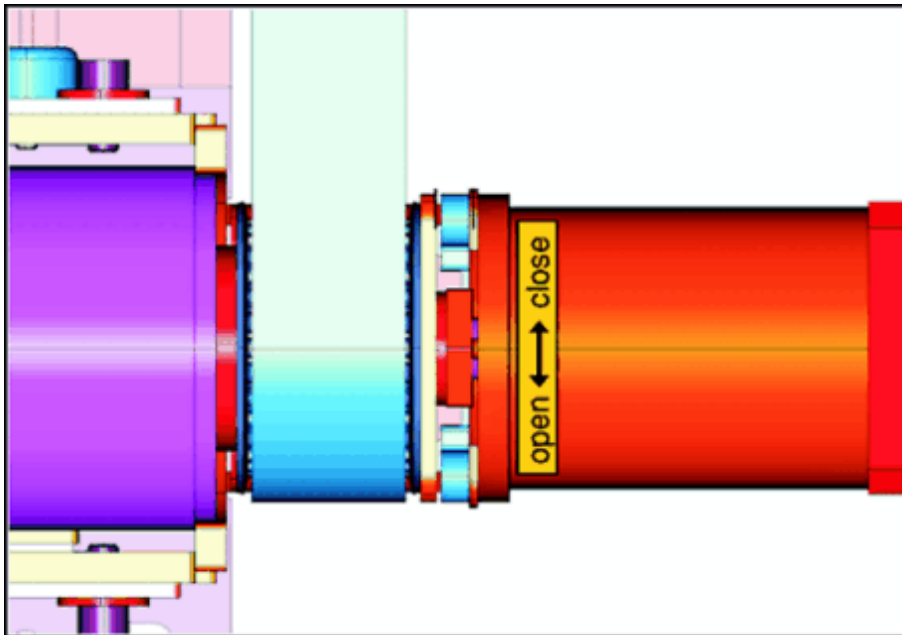



Fig. 185: *M6 hex screws of MDSD bracket*

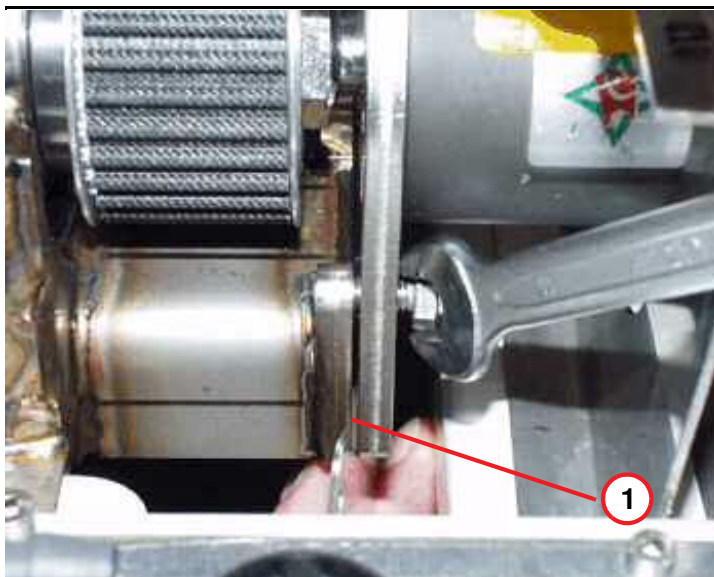

- While torquing the 22 mm nut, a small clearance appears.

(1) Shim gap appears when securing the 22 mm nut

Fig. 186: Insert shim plates

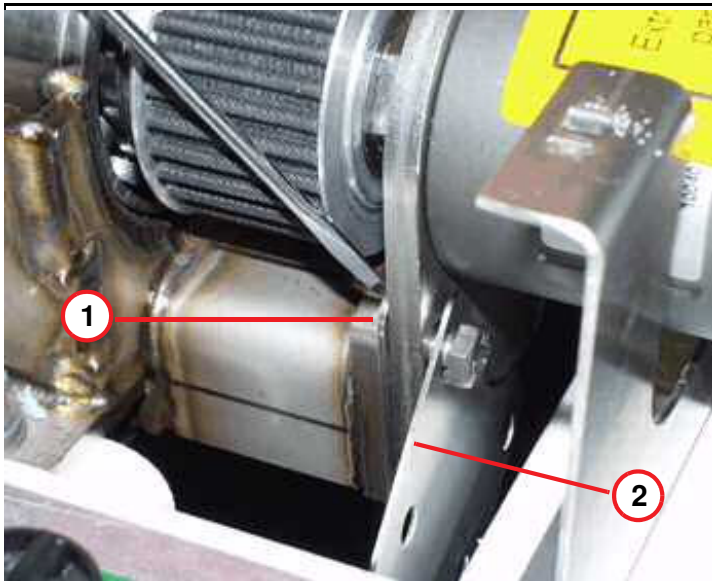

- (1) Shim plates to be added
(2) Spare shim plates

- Use the shim plates to fill the clearance between the motor carrier and the axle (delivered with the spare unit, or component of the original part). Tighten all four hex screws.

3.8.5.2 MDSD drive of partial drive, not valid for Trio A Tim Systems




CAUTION

Strong magnetic field!

Magnetic objects may cause injury.

- » Hold the ferromagnetic part with both hands while moving it to the correct position and to remove the unit directly out of the field!


Fig. 187: MDSD clutch removed

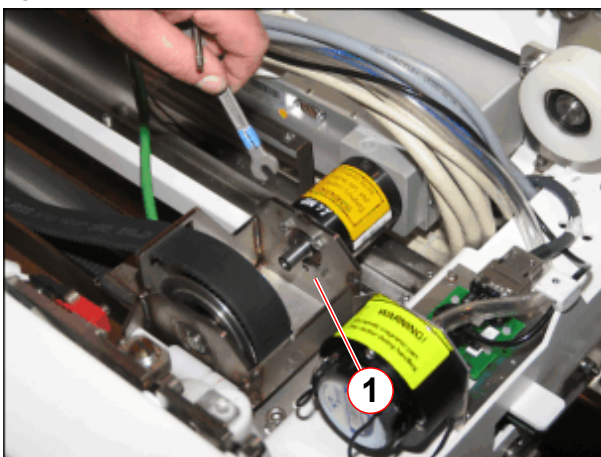

- (1) Mounting screws of MDSD

- Remove the hex screws (1) of the mounting bracket from the MDSD.


Fig. 188: Clutch of MDSD, partial drive

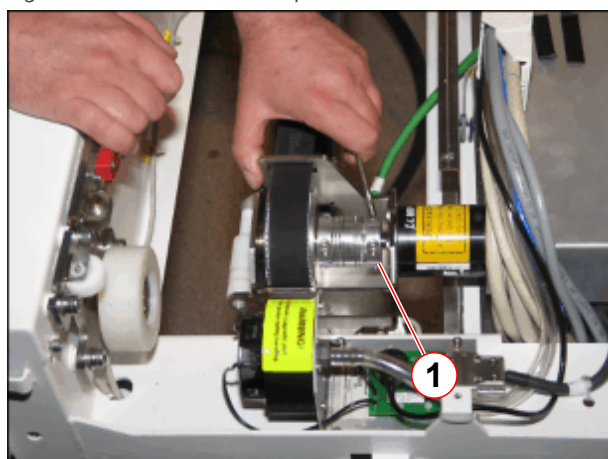

(1) Allen screws of MDSD clutch

- Remove Allen screws of MDSD clutch (1).
- Hold the ferromagnetic part with both hands while moving it to the correct position and to remove the unit directly out of the field!


Fig. 189: Assembling the partial drive MDSD

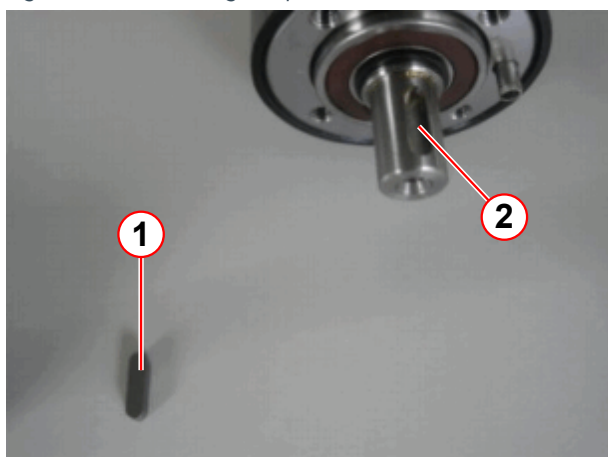

(1) Key
(2) Notch of the MDSD axis


- For partial drive tables insert the key (1), delivered with the spare part, into the notch of the MDSD axis.

Fig. 190: Key inserted in MDSD axis

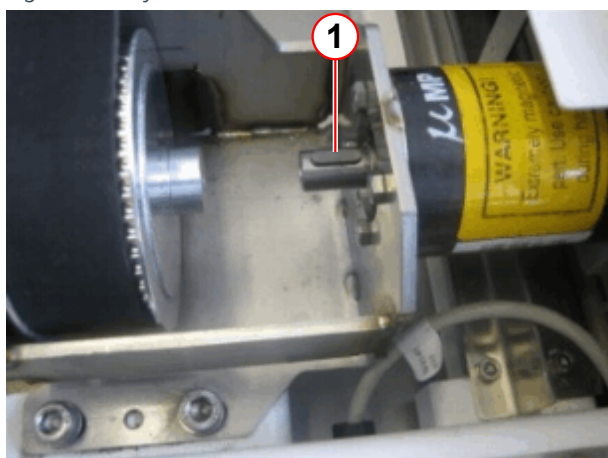

(1) Key


Fig. 191: Fixation of the clutch

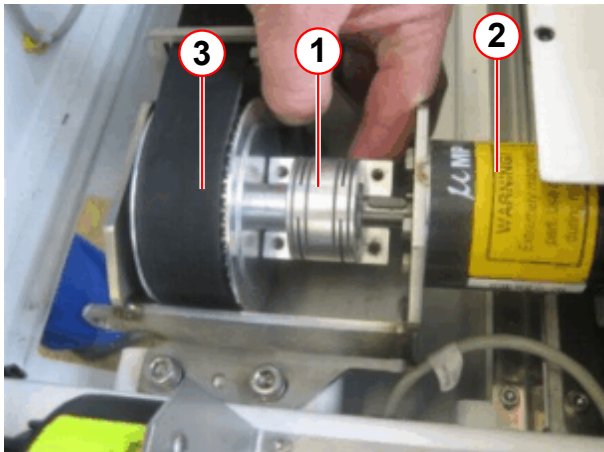

- (1) Clutch
- (2) MDSD for partial drive
- (3) Drive unit TK


- Place the clutch as shown in pos. (1).

Fig. 192: Fixation of the clutch

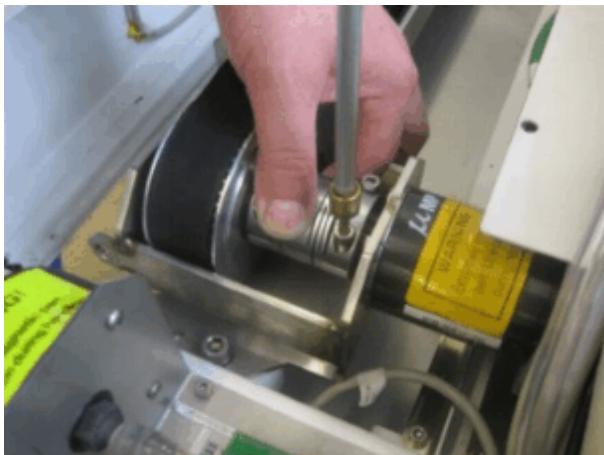


- Position the other half of the clutch so that the key fits into the notch of the clutch and tighten the screws of the clutch.

Fig. 193: Key position with clutch

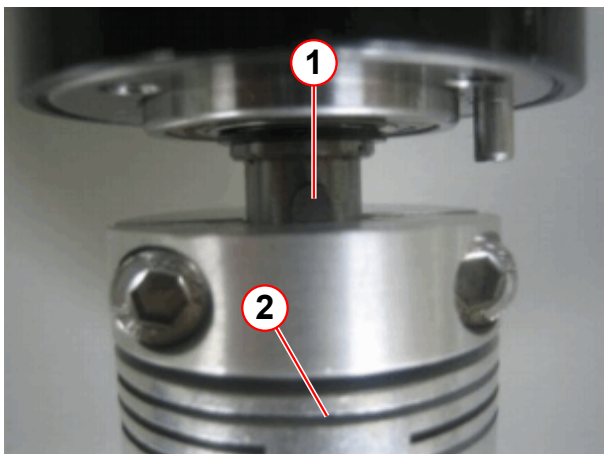

- (1) Key
- (2) Clutch

Fig. 194: Clutch installed

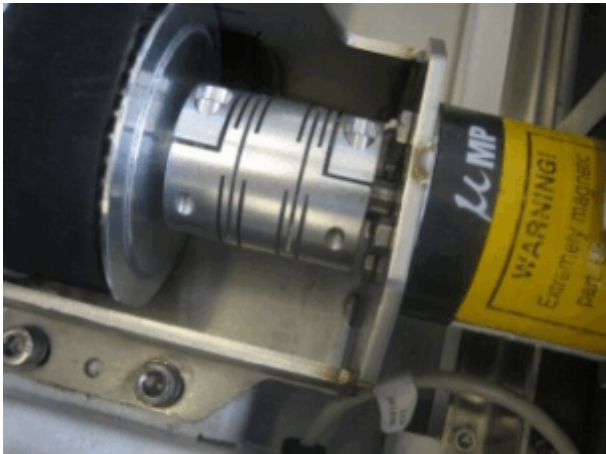


For further assembling of the partial drive, please refer to (→ Assembly / Page 167), not valid for Trio A Tim Systems.

3.8.6 Adjustments

For adjustment of the toothed belt, please refer to (→ Adjustments / Page 170), not valid for Trio A Tim Systems.

3.8.7 Concluding steps

Fig. 195: MDSD FOC connections

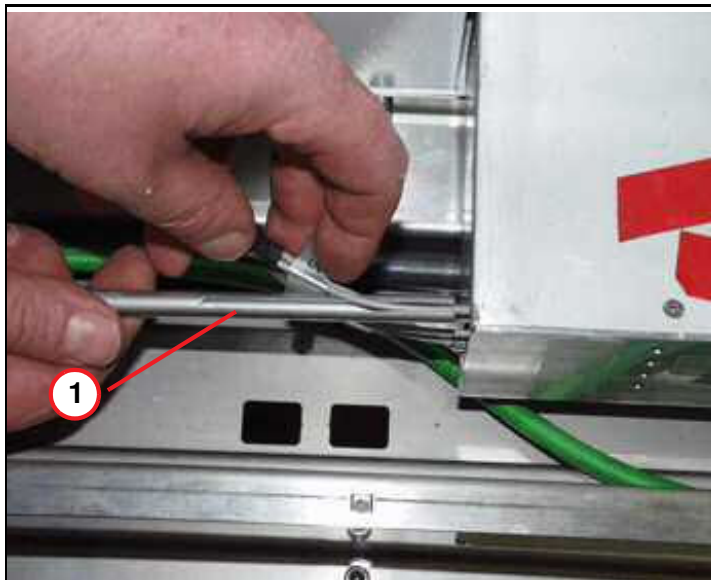

(1) FOC feed-in tool

- Reconnect the cables and FOCs on the MDSD drive and use the feed-in tool (service aid) delivered with the spare part.

- Reinstall the covers in reverse order, tighten all screws at the end of the assembling procedure.

- Switch ON F22 (RF-room).

3.8.8 Adjustments

For adjustment of the toothed belt, please refer to (→ Adjustments / Page 170), not valid for Trio A Tim Systems.

3.8.9 Final checks

- Please refer to final checks horizontal drive unit (→ Final checks / Page 157).

3.9 Horizontal drive unit

3.9.1 General

The horizontal drive unit is responsible for the horizontal movement of the ptab. The unit can only be replaced in full.

There are two types of horizontal drives unit in the patient table.

1. "Telescopic Drive" that translates a relatively small movement range into a wide range to support whole body applications.
2. "Partial drive" allows the horizontal table movement without whole body application and a limited movement range (not valid for Trio A Tim Systems).

3.9.2 Tools and time required

3.9.2.1 Time

Approx. 120 min

3.9.2.2 Tools

Non-magnetic tool kit

Philips screwdriver

3.9.3 Preliminary

- Switch OFF F22 (RF-room) in the ACC



According to the system/room specifications the movement range of the horizontal movement can be reduced.

- If the movement range of the horizontal movement is reduced it may be necessary to remove the end stops. (→ Fig. 289 Page 173)

3.9.4 Disassembly

- » Please follow the procedure if the horizontal drive unit needs to be replaced.

**CAUTION**

The horizontal drive units consist of magnetic objects!

Magnetic objects may cause injury.

- » Maintain a safe distance to the magnetic field!

3.9.4.1 Tabletop plate



Be careful and avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.



The tabletop plate may fall - support the tabletop plate.

Basic table

Fig. 196: Coil plug cover (foot end)



- Remove the Allen screws on top of the cover from the coil plugs.


Fig. 197: Remove Allen screws



- Remove the lower Allen screws of the cover from the coil plugs.

Removable table


Fig. 198: RF-Cover RTT

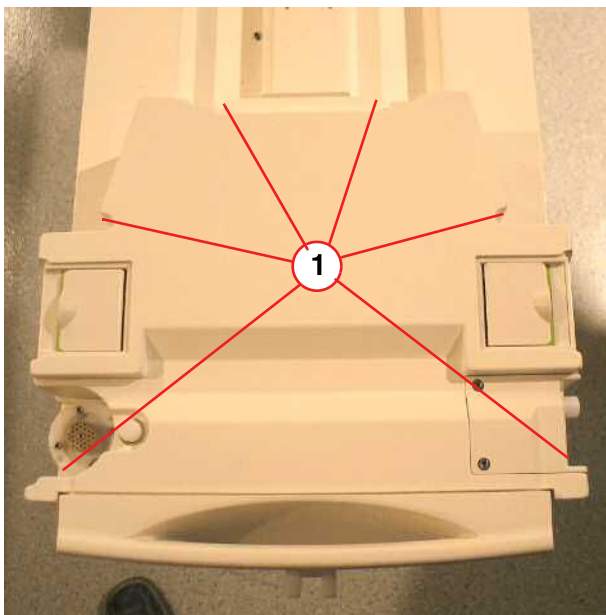

(1) Screws to be removed

- Remove the carrier plate of the RTT.
- Remove six screws on top of the cover from the coil plugs.
- Disconnect cables.

Fig. 199: Additional mounting bracket for foot end cover (newer systems)

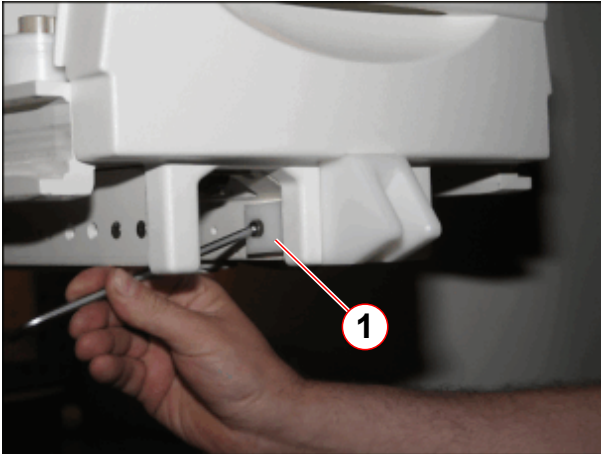

(1) Allen screw


- On newer systems, remove Allen screw of mounting bracket.

Fig. 200: Foot-end cover



- Remove screws of the foot-end cover.


Fig. 201: Foot end cover, removed

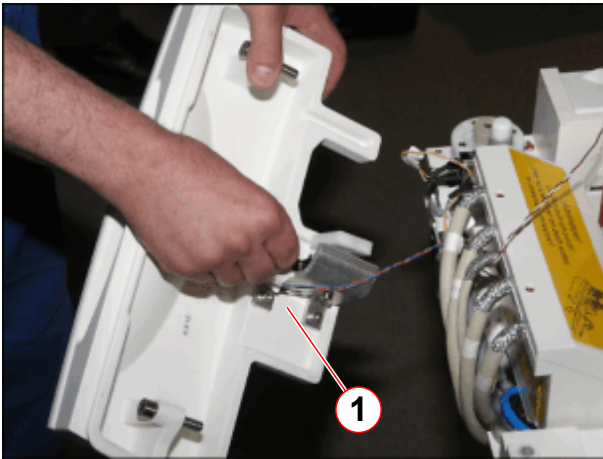

(1) Disconnect cables

- Disconnect cable connections and remove the foot-end cover.
- Slide out the lower cover on the bottom of the patient table.

Fig. 202: Mechanical end stop of partial drive unit

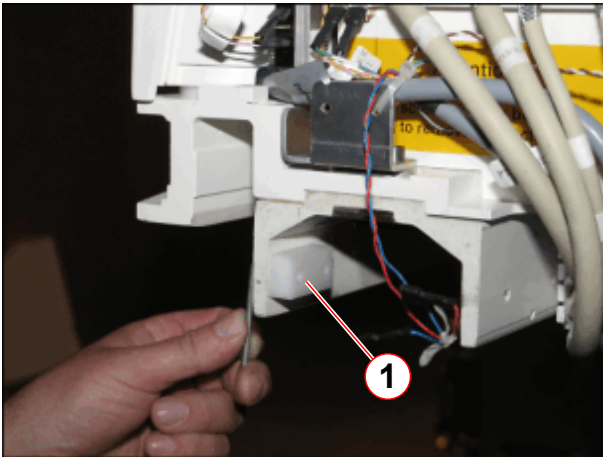

(1) End stop to be removed

- Remove mechanical end stop.

Fig. 203: Energy chain




- Remove the lower part of the energy chain

Whole Body / telescopic drive


- For the telescopic drive: Release the emergency clutch and carefully push the tabletop plate manually into the magnet until the cantilever arm is no longer touching the tabletop plate. Support the table to ensure it does not fall.

Fig. 204: Tabletop, magnet end

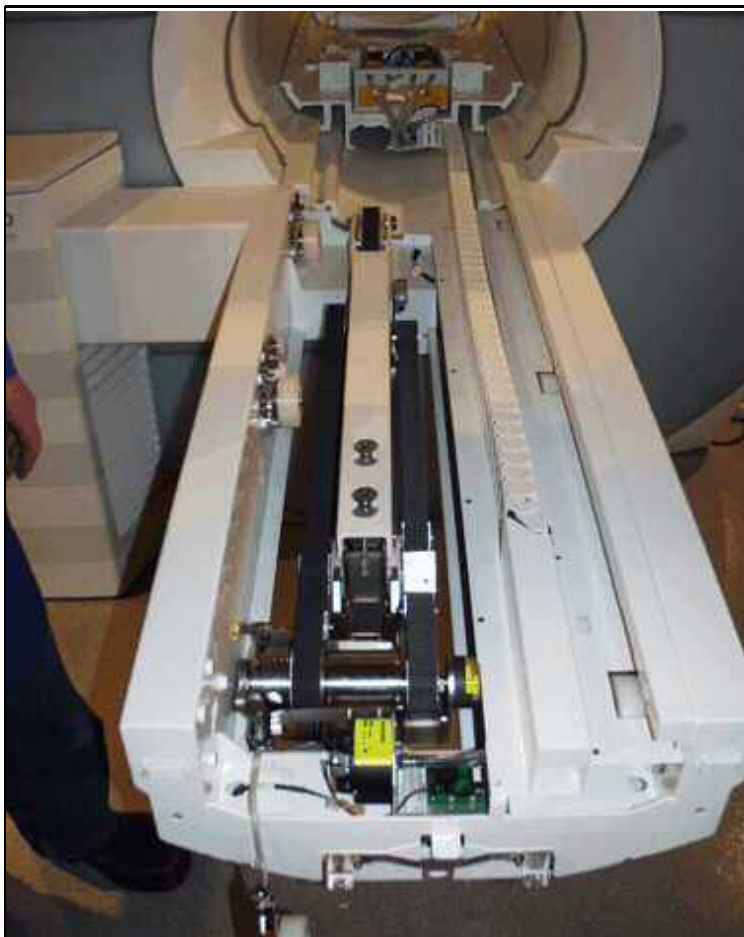

Partial drive, not valid for Trio A Tim Systems


Fig. 205: Tension unit for toothed belt

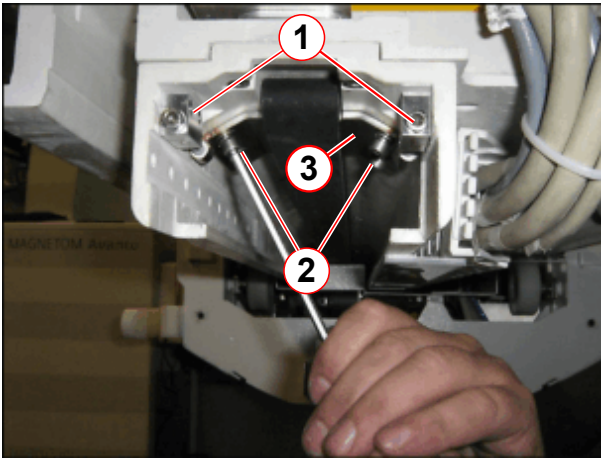

- (1) Allen screws for adjusting belt tension
- (2) Attachment screws
- (3) clamping plate

- Release tension of toothed belt (1), remove Allen screws (2) and tooth belt tension unit.


Fig. 206: Cover of toothed belt mount (head end)

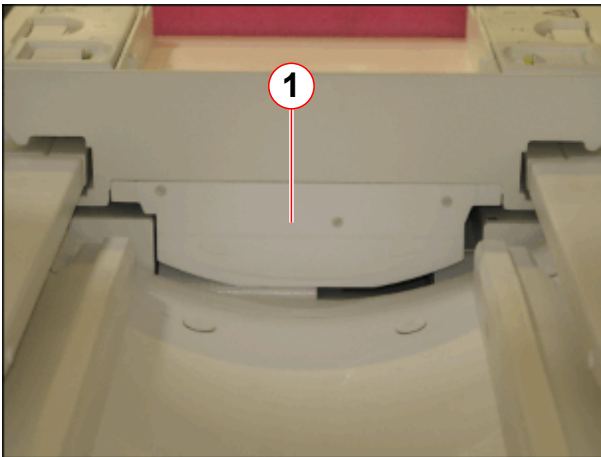

- (1) Cover of toothed belt mount

- Remove the cover of the toothed belt mount of toothed belt (1) on the head end from the ptab.


Fig. 207: Toothed belt attachment (head end)

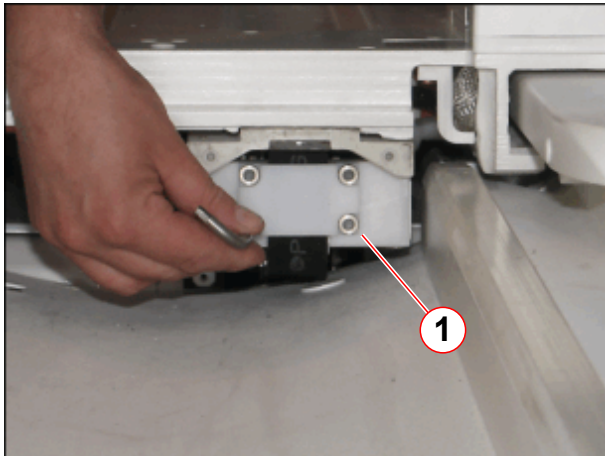

(1) Allen screws of toothed belt mount

- Remove the toothed belt mount (1) on the head end from the ptab.

Fig. 208: Tabletop moved into the magnet

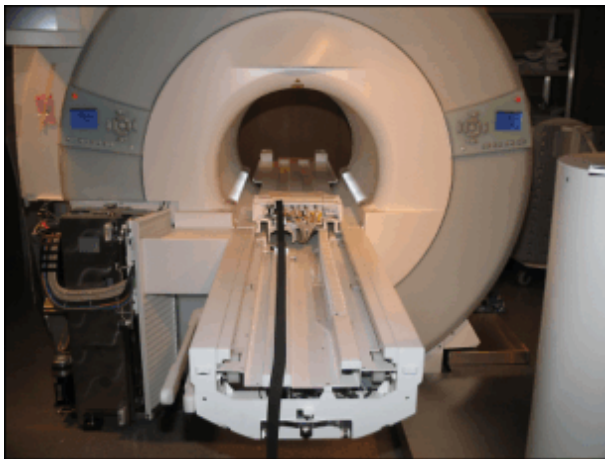


- Carefully push the tabletop plate manually into the magnet.

3.9.4.2 Covers / cables / energy chain



Be careful and avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.


Fig. 209: Plastic inner cover

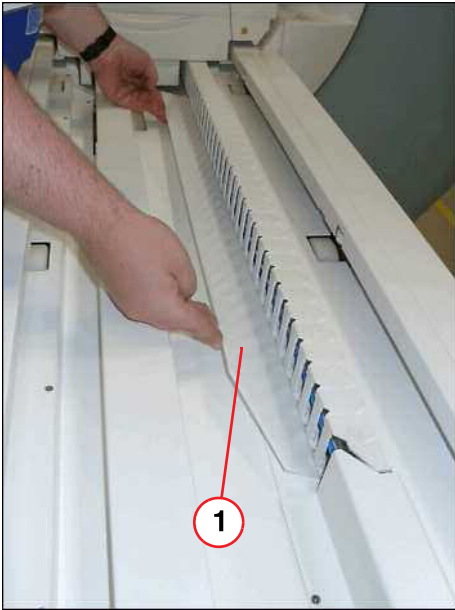

(1) Plastic inner cover

- Remove inner plastic cover.

Fig. 210: Remove cover of energy chain

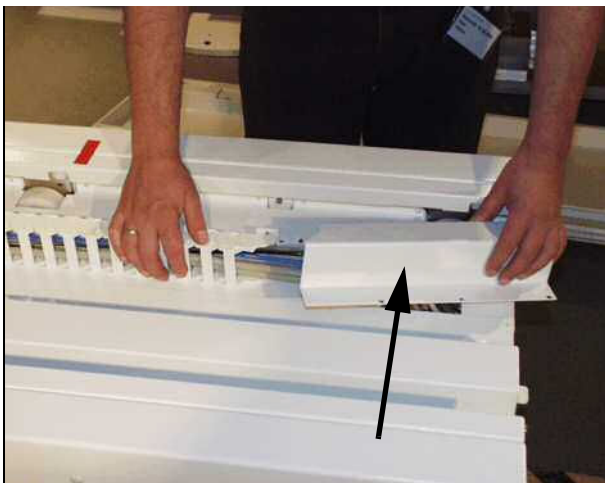


- Remove the upper energy chain cover.

Fig. 211: Remove energy chain



- Remove the upper energy chain of the bracket.


Fig. 212: Inner main cover

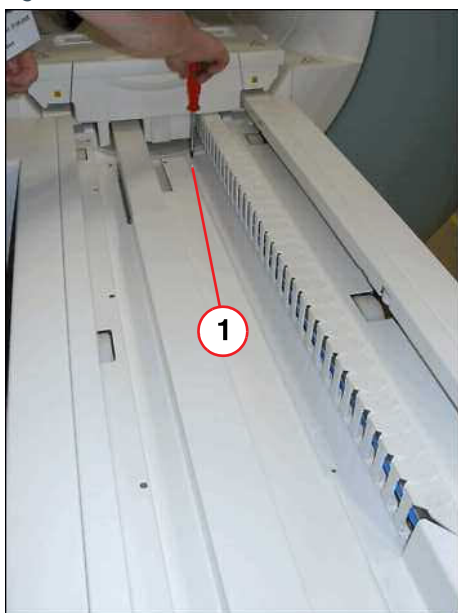

(1) Screws to be removed

- Remove screws of main inner cover.

Fig. 213: Allen screws

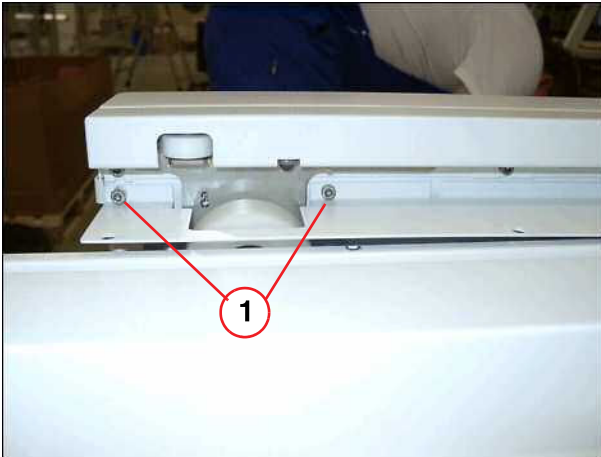


(1) Screws to be removed


- Remove two Allen screws at the left side of the table end (foot) before you remove the main inner cover.

- Remove the inner main cover.

Fig. 214: Inner main cover

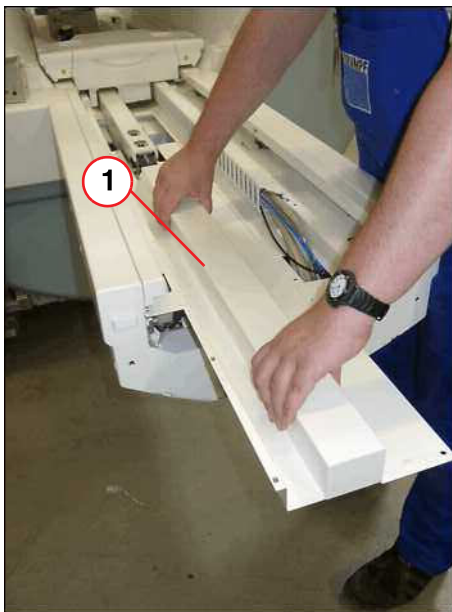


(1) Inner main cover

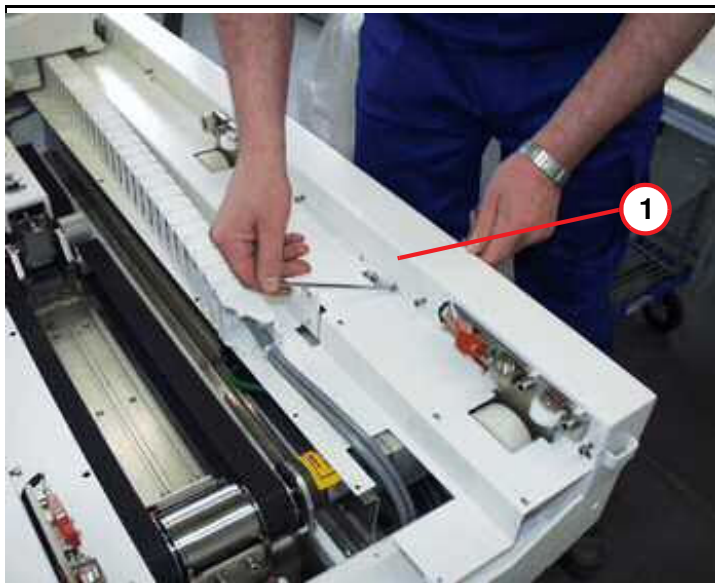Fig. 215: Inner cover

(1) Top inner cover


- Remove the right guide rail cover and the inner top cover.

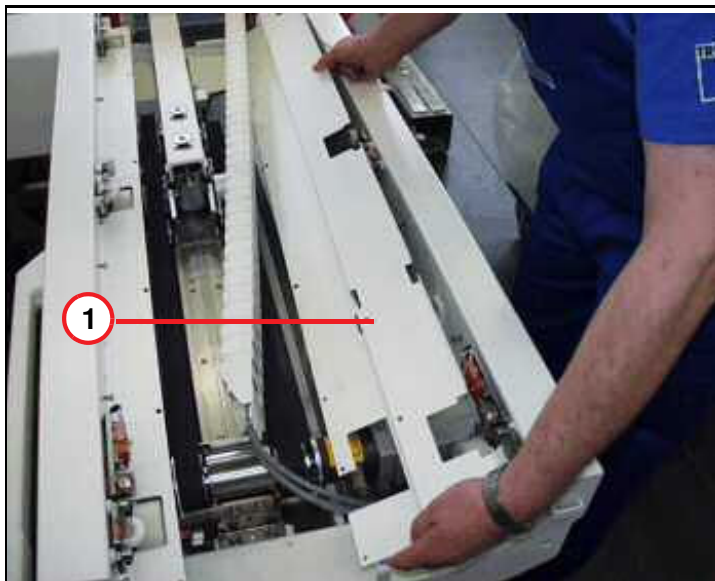Fig. 216: Inner cover

(1) Inner cover, right side

- Remove the inner right cover.


Fig. 217: Inner cover, left side

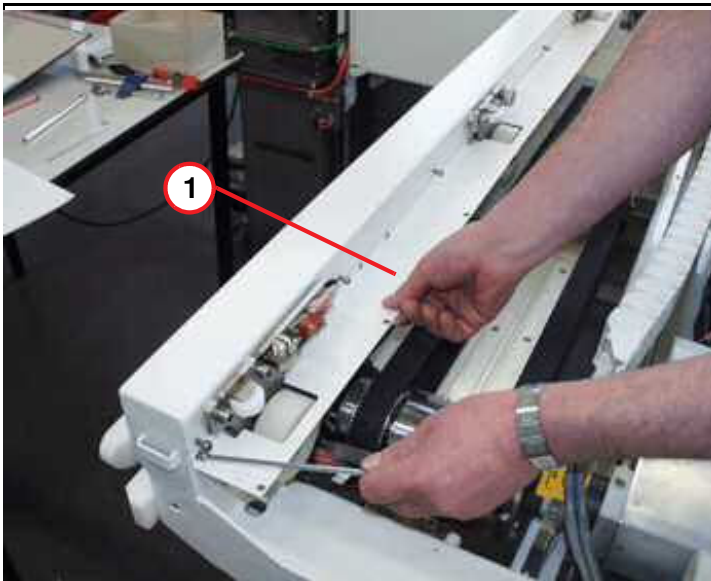

(1) Inner cover, left side

- Remove the left guide rail and the inner cover, left side.

Fig. 218: Covers removed

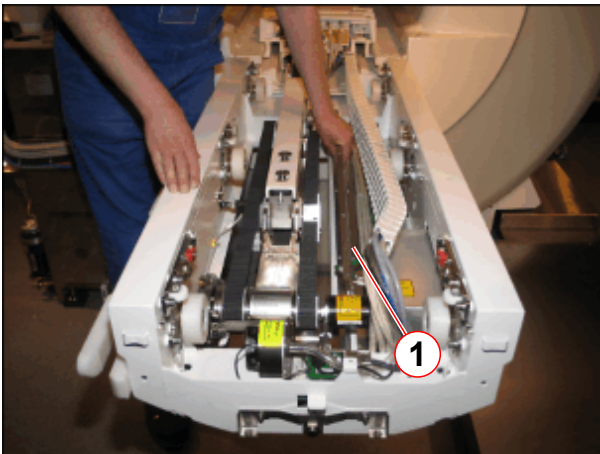

(1) Support plate

- Remove support plate.

3.9.4.3 Removal of horizontal drive unit / telescopic drive




CAUTION

The horizontal drive units consist of magnetic objects!

Magnetic objects may cause injury.

- » Maintain a safe distance from the magnetic field.
- » Fasten the horizontal drive unit with the arresting cable.


Fig. 219: Control switches

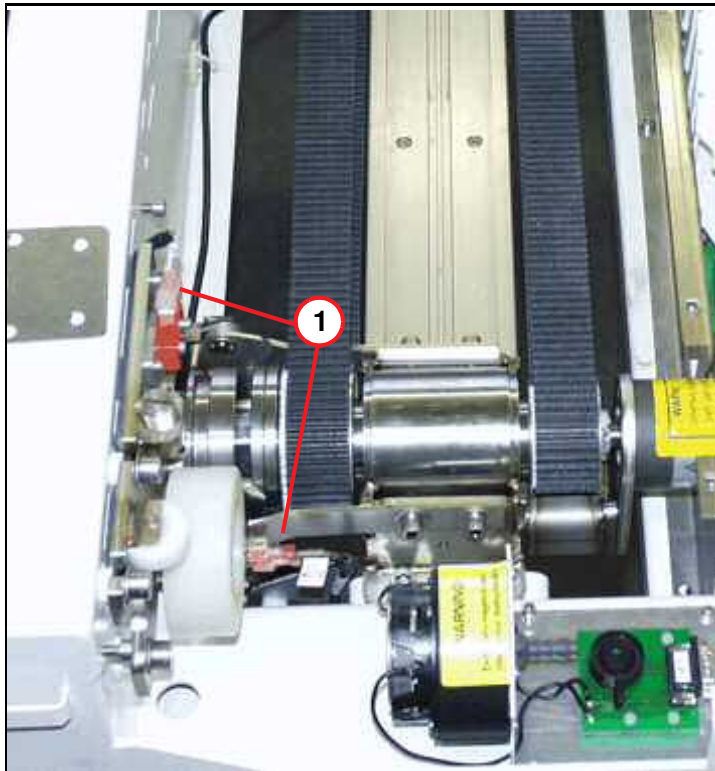

(1) Switch S51 and S32

- Disconnect the switches, S51 (emergency release), S32.

Fig. 220: Plug S10

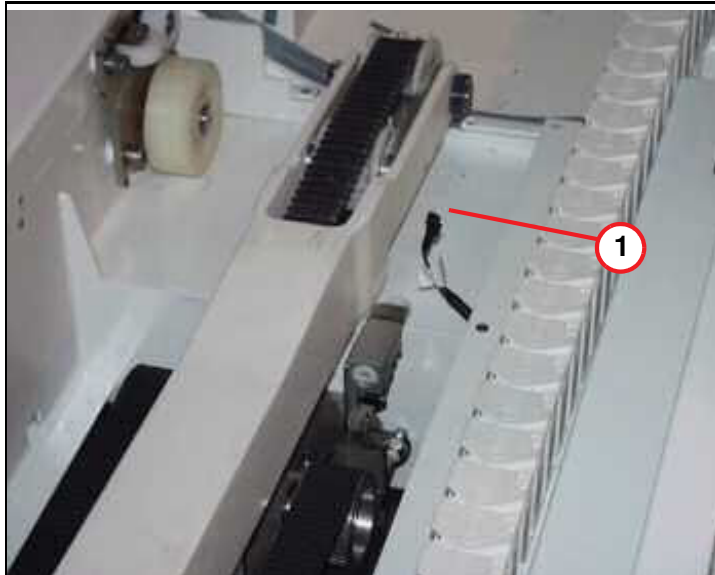


(1) S10 plug, disconnected

- Disconnect the cables of S10 (reference switch).


- For removing and mounting, the horizontal drive unit (p/n 8112257) must be secured with the arresting cable.
- Secure the horizontal drive with the arresting cable as shown below

Fig. 221: Arresting cable secured on the horizontal drive

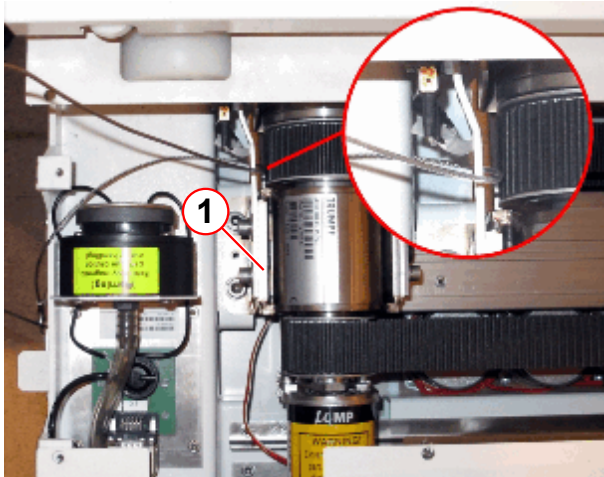

(1) Horizontal drive bracket

Fig. 222: Arresting cable secured on the mounting angle

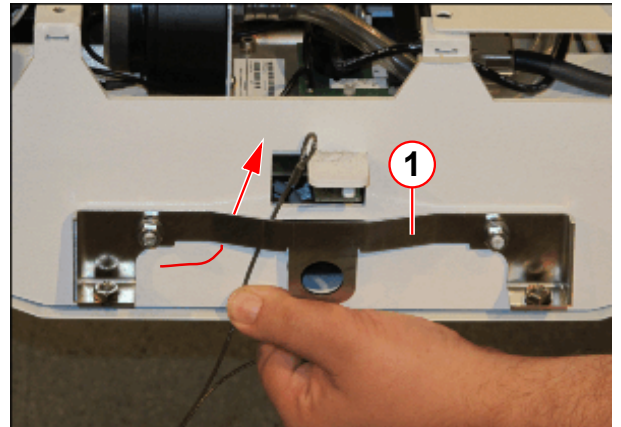

(1) Mounting angle


Fig. 223: Hook



(1) Loops
(2) Hook

- Secure the arresting cable loops (pos.1) by closing the hook (pos.2).

Fig. 224: Position of drive bracket

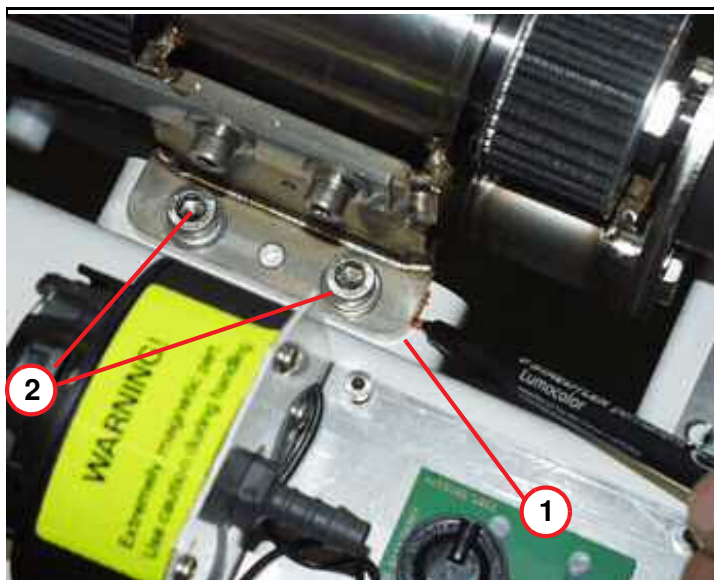

(1) Original position of drive bracket
(2) Fastener of horiz. drive unit


- Mark the position of the carrier plate with respect to the drive bracket.
- Remove the two M8 Allen screws of the horizontal drive bracket.

Fig. 225: M8 Allen screw, horizontal drive unit



- Remove the M8 Allen screw underneath the patient table.


Fig. 226: Example - removing FOCs from MDSD

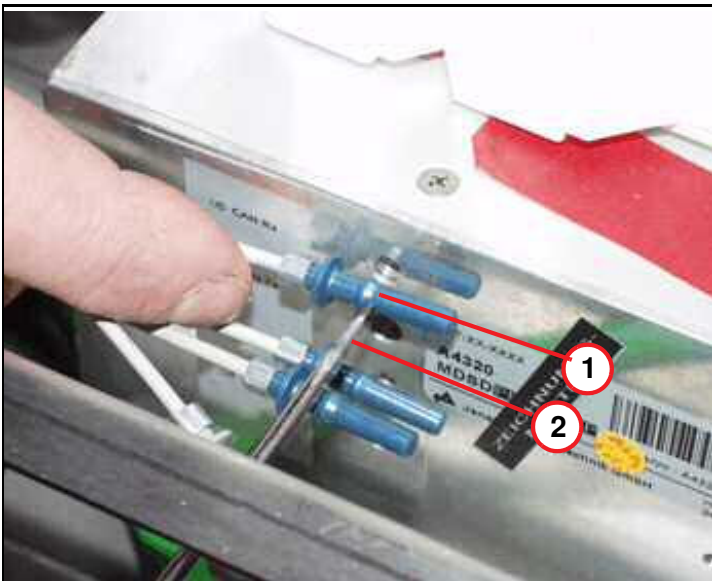

- (1) Cam of FOC connector
- (2) Screwdriver, flat blade

- The figure below shows the FOC plug connection (principle) inside the MDSD device. Use a flat blade screwdriver to disconnect the FOCs inside the MDSD device.

Fig. 227: Fiber-optic cables of MDSD

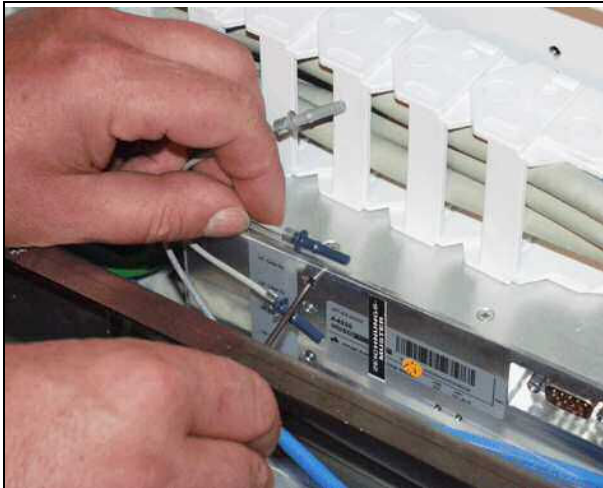

- Feed the flat blade screwdriver through the hole in the aluminum cover and into the FOC's plug connection. Try to reach the cam of the FOC and carefully push it sideways. At the same time, pull it out from the FOC. Refer to the illustration.



CAUTION

Due to the strong magnetic parts of the horizontal drive, it is very important to move the horizontal drive always away from the magnet!

Do not move the drive closer to the magnet than absolutely necessary.

- » Maintain a safe distance from the magnetic field.


Fig. 228: Horizontal drive removed on the service platform



- Store the secured horizontal drive at the service platform at the patient end of the table.
- Replace the horizontal drive with the new one. Replace the original MDSD device as described in the MDSD chapter. (→ MDSD drive of whole body drive / Page 116)
- For the replacement, the arresting cable must be used.

3.9.4.4 Removal of partial drive unit, not valid for Trio A Tim Systems


Fig. 229: Partial drive assembly

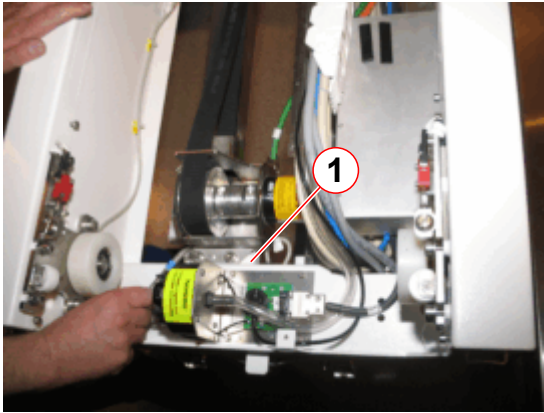

(1) Allen screws of mounting bracket

- Remove Allen screws of partial drive unit, foot end side.


Fig. 230: Partial drive unit, button view

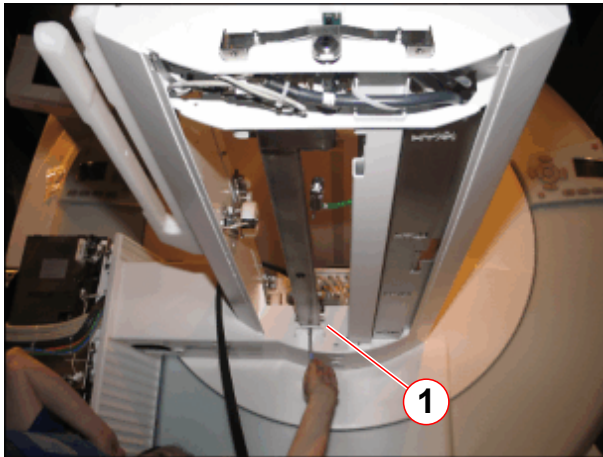

(1) Allen screw for fastening the partial drive unit

- Remove Allen screw on the bottom of the partial drive unit.


Fig. 231: Example - removing FOCs from MDSD

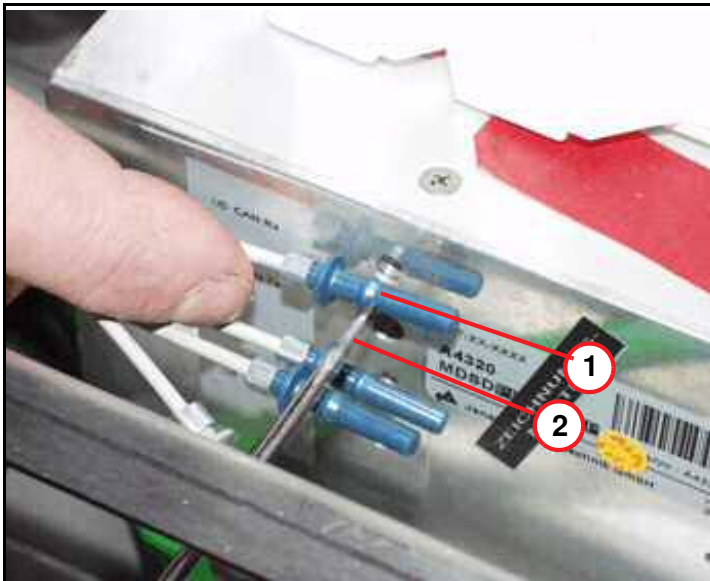

- (1) Cam of FOC connector
- (2) Screwdriver, flat blade

- The figure below shows the FOC plug connection (principle) inside the MDSD device. Use a flat blade screwdriver to disconnect the FOCs inside the MDSD device.

Fig. 232: Fiber-optic cables of MDSD

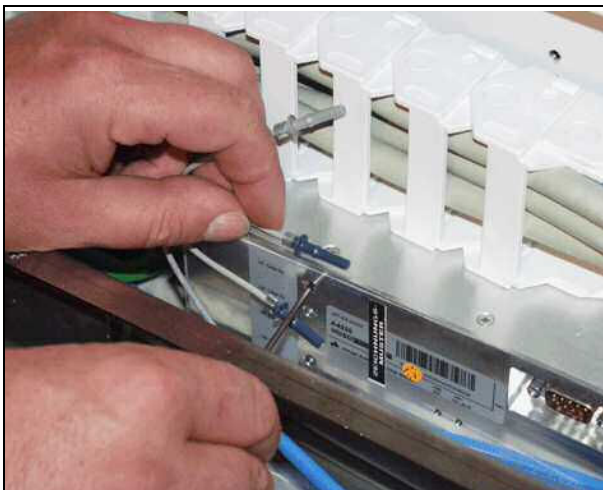

- Feed the flat blade screwdriver through the hole in the aluminum cover and into the FOC's plug connection. Try to reach the cam of the FOC and carefully push it sideways while at the same time pulling it out from the FOC. Refer to illustration.

Fig. 233: Partial drive assembly, removed



(1) Position of toothed belt

- Remove the partial drive unit and store it on a clean workbench.
- Replace the original MDSD device as described in the MDSD chapter. (→ MDSD drive of whole body drive / Page 116)

3.9.5 Assembly

3.9.5.1 Horizontal drive unit / telescopic drive



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.



CAUTION

Due to the strong magnetic parts of the horizontal drive, it is very important to move the horizontal drive always away from the magnet!

Do not move the drive closer to the magnet than absolutely necessary.


» **Maintain a safe distance from the magnetic field.**

Fig. 234: Horizontal drive removed on the service platform




- Store the secured horizontal drive at the service platform at the patient end of the table.
- Secure the horizontal drive with the arresting cable as shown in the pictures below.
- Replace the horizontal drive with the new one. Replace the original MDSD device as described in the MDSD chapter. (→ MDSD drive of whole body drive / Page 116)

Fig. 235: Arresting cable secured on the horizontal drive

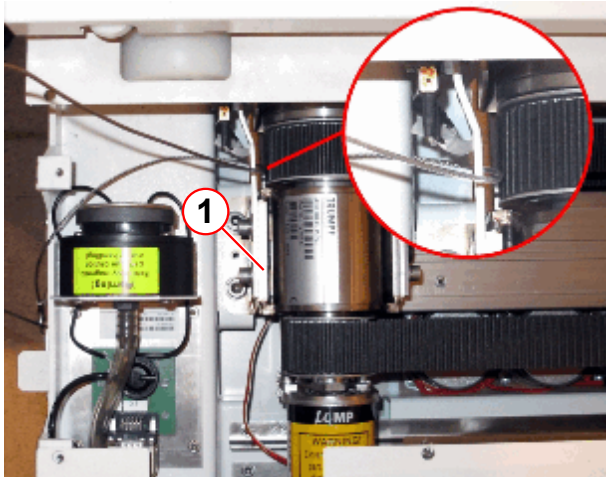


(1) Horizontal drive bracket

Fig. 237: Hook

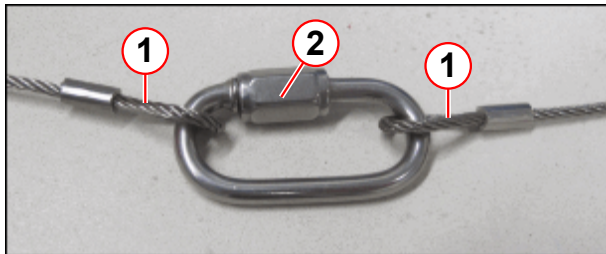


(1) Loops
(2) Hook

Fig. 238: Assembling the partial drive MDSD

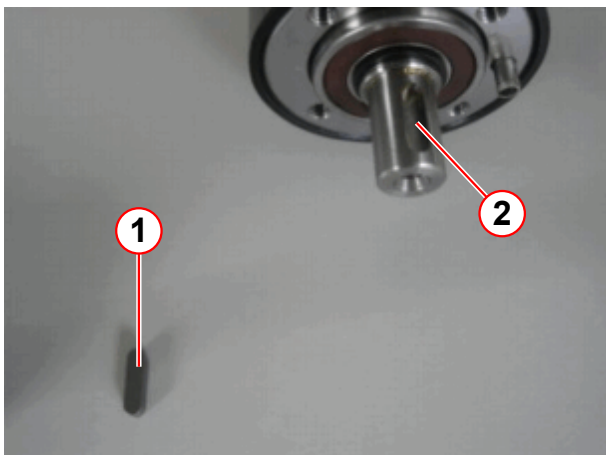

(1) Key
(2) Notch of the MDSD axis


Fig. 236: Arresting cable secured on the mounting angle

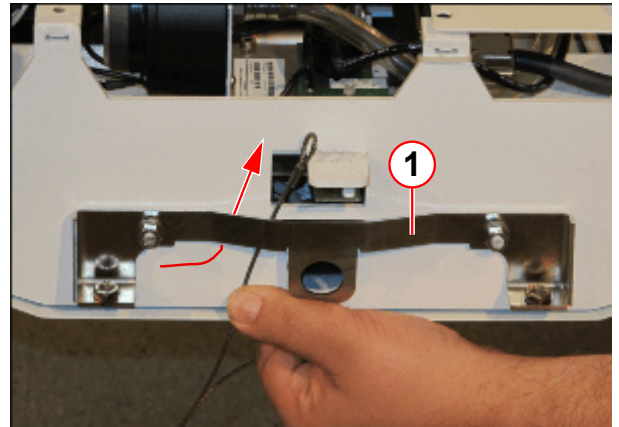

(1) Mounting angle


- Secure the arresting cable loops (pos.1) by closing the hook (pos.2).
- For partial drive tables insert the key (1), delivered with the spare part, into the notch of the MDSD axis.

Fig. 239: Key inserted in MDSD axis

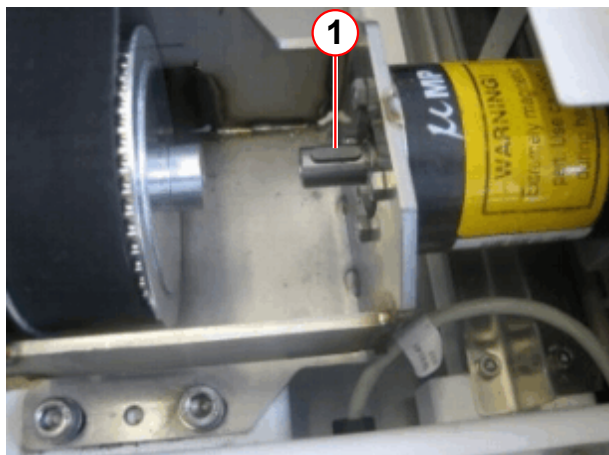

(1) Key


Fig. 240: Fixation of the clutch

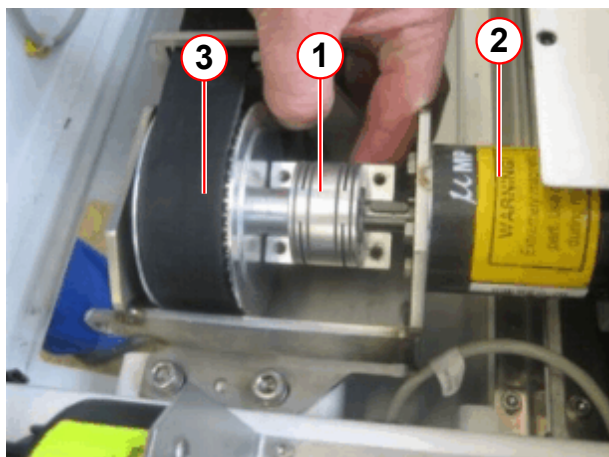

(1) Clutch
(2) MDSD for partial drive
(3) Drive unit TK


- Place the clutch as shown in pos. (1).

Fig. 241: Fixation of the clutch

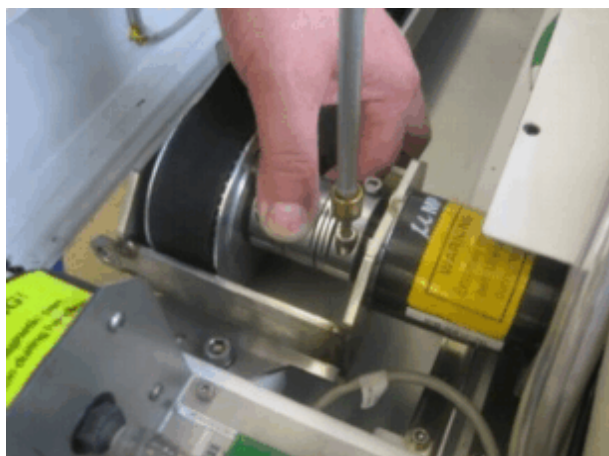

- Position the other half of the clutch so that the key fits into the notch of the clutch and tighten the screws of the clutch.

Fig. 242: Key position with clutch

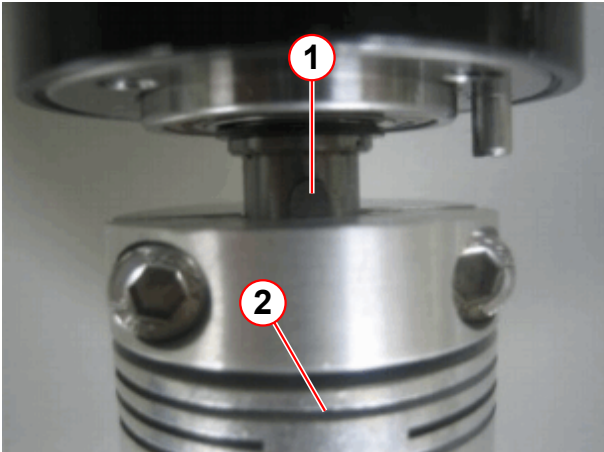

- (1) Key
- (2) Clutch

Fig. 243: Clutch installed


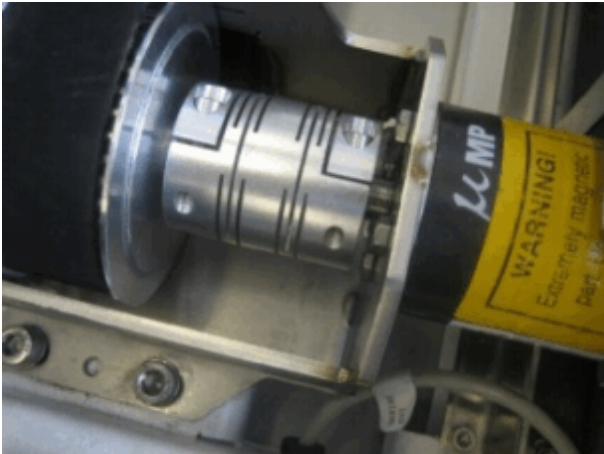


Fig. 244: Position of drive bracket

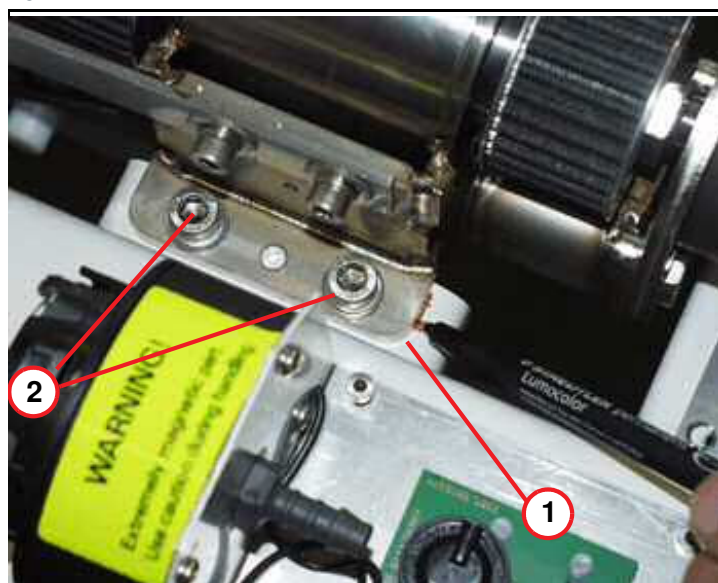

- (1) Original position of drive bracket
(2) Fastener of horiz. drive unit

- Replace the horizontal drive with the new one and use the key delivered with the spare part as described before. Replace the original MDSD device as described in the MDSD chapter. (→ MDSD drive of whole body drive / Page 116)
- Reinstall the horizontal drive unit and adjust the drive bracket to the marked position.
- Reinsert the Allen screws of the horizontal drive mechanism, apply Loctite 221.


Fig. 245: M8 Allen screw, horizontal drive unit



- Reinsert the Allen screw of the horizontal drive underneath the table, apply Loctite 221.


- Removal the arresting cable.
- Reconnect switch S10 and S51.
- Reconnect all cables to the MDSD device, use the feed-in tool as described in the MDSD chapter.
- Reinstall the assembly of the guide roller, apply Loctite 221.
- Reconnect switch S32.
- Reinstall all inner table covers.

Fig. 246: Distance of cantilever arm

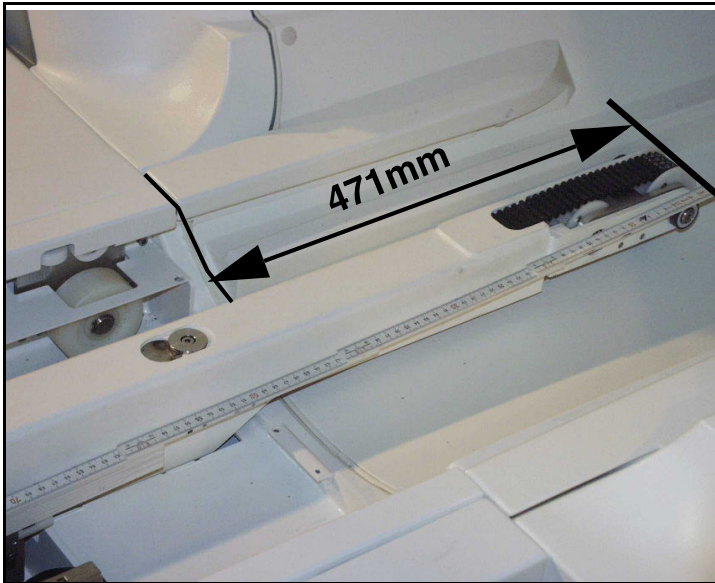

- Before you insert the new tabletop plate, please do a check of the distance from the end of the cantilever arm to the support frame. Adjust to 471 mm.
- Carefully connect the tabletop plate to the toothed belt of the horizontal drive.

Fig. 247: Mark on the telescopic arm

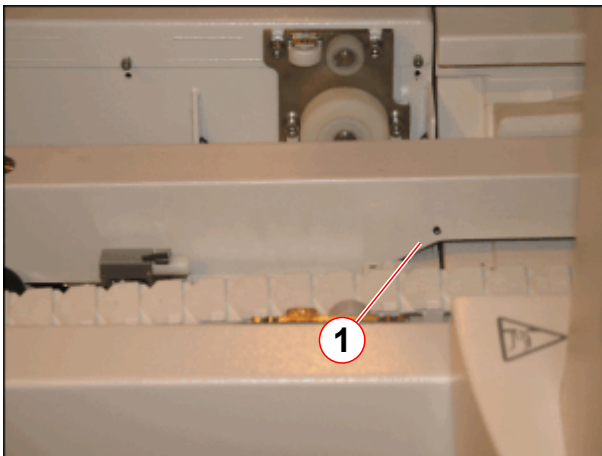

(1) Marking hole in the telescopic arm

- On newer systems the distance mark shows the correct length of the cantilever arm.

3.9.5.2 Partial drive unit, not valid for Trio A Tim Systems



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.


Fig. 248: Partial drive assembly, removed



(1) Position of toothed belt

- Reinsert the horizontal drive unit.


Fig. 249: Tabletop into the magnet; toothed belt route

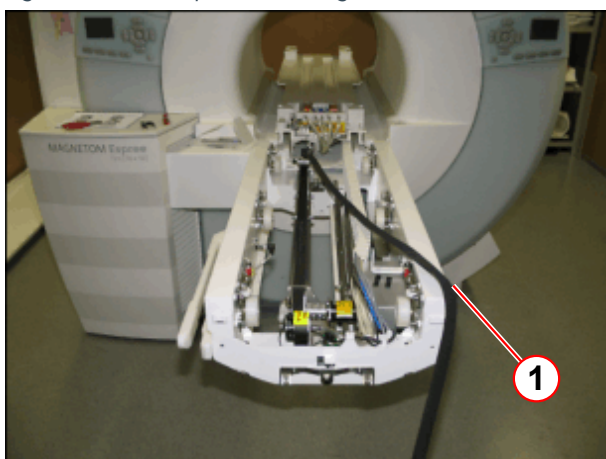

(1) To avoid misassembling, lay toothed belt over side

- Before you insert the the cover, ensure the correct toothed belt route.


Fig. 250: Partial drive assembly

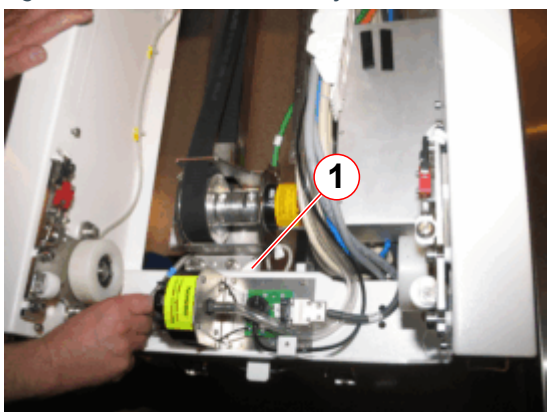

(1) Allen screws of mounting bracket

- Reinstall the drive unit.
- Reinsert Allen screws of partial drive unit, foot end side, apply Loctite 221.

Fig. 251: Partial drive unit, button view

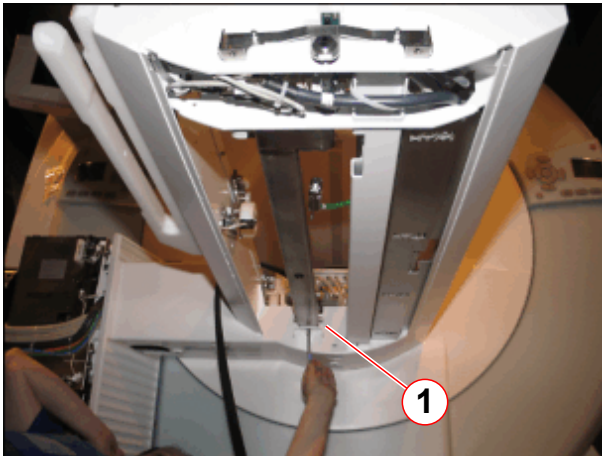


(1) Allen screw for fastening the partial drive unit

- Reinsert Allen screw on the bottom of the partial drive unit, apply Loctite 221.


- Reconnect switches.
- Reconnect all cables to the MDSD device, use the feed-in tool as described in the MDSD chapter. (→ MDSD drive of whole body drive / Page 116)
- Reinstall all inner table covers.

Fig. 252: Tabletop moved into the magnet

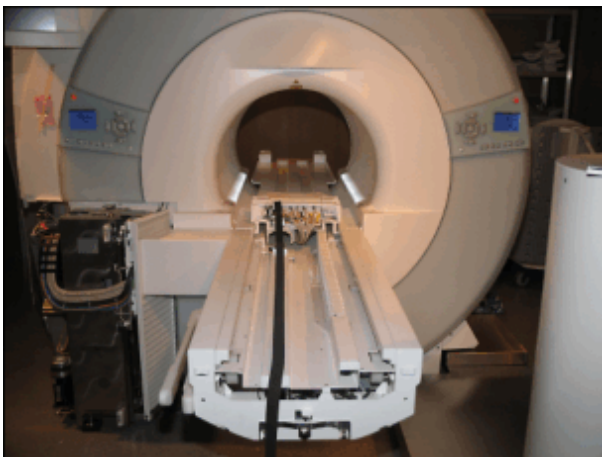

- Make sure that the toothed belt runs correct under the table top.


Fig. 253: Toothed belt attachment (head end)

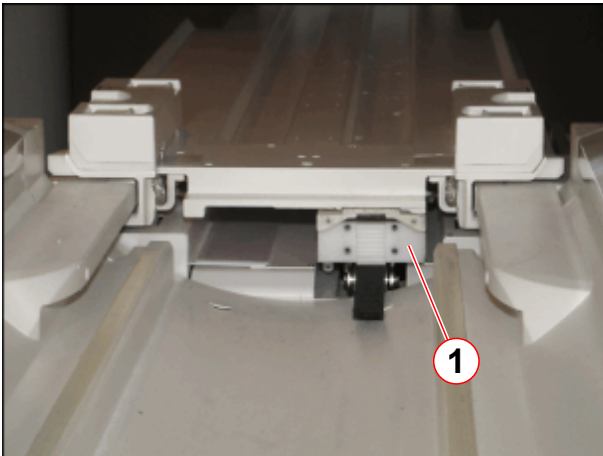

(1) Toothed belt mount (head end)

- Carefully pull the tabletop plate manually back to the food end side (patient end).
- Install the mount of toothed belt on the head end from the ptab.

Fig. 254: Cover of toothed belt mount (head end)

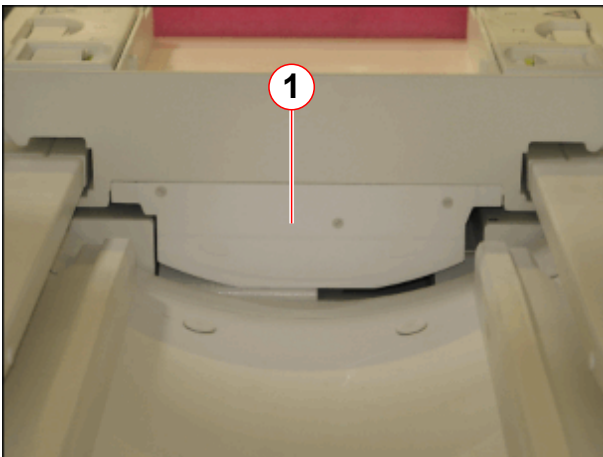

(1) Cover of toothed belt mount


- Install the cover of the toothed belt mount(1) on the head end from the ptab.

Fig. 255: Mounting the toothed belt (foot end)



- Principle of toothed belt tension unit, located on the foot end side of the ptab.

Fig. 256: Service tool for toothed belt tension

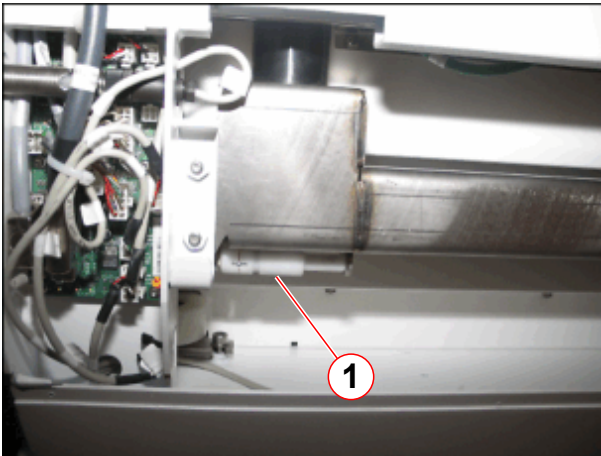

(1) Storage position of service tool for toothed belt tension

- The service aid for adjustment of the toothed belt tension is located underneath the table top.
- Snap out the service aid for adjustment of the toothed belt tension.

3.9.6 Concluding steps

- Install the energy chain.



Use the same length of screws when reinstalling the energy chain


- Make sure you replace the covers in the proper sequence and tighten all screws at the end of the assembling procedure.
- Install the end stops.
- Install the remaining covers.
- Install foot-end plate and cover from the coil plug.

3.9.7 Adjustments, not valid for Trio A Tim Systems

Fig. 257: Service tool, toothed belt tension




Fig. 258: Tension unit for toothed belt

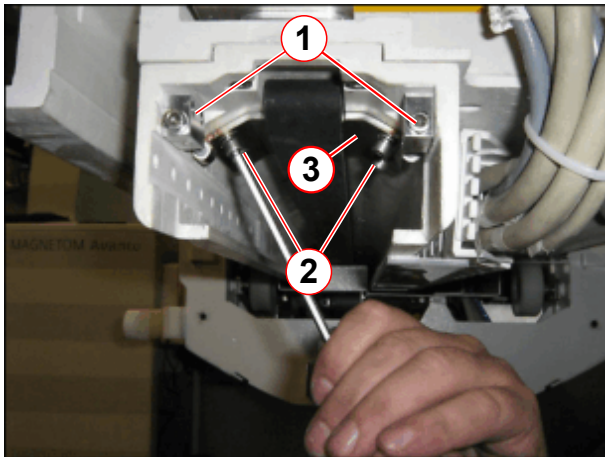

- (1) Allen screws for adjusting belt tension
- (2) Attachment screws
- (3) clamping plate

Fig. 259: Service tool for toothed belt tension in use

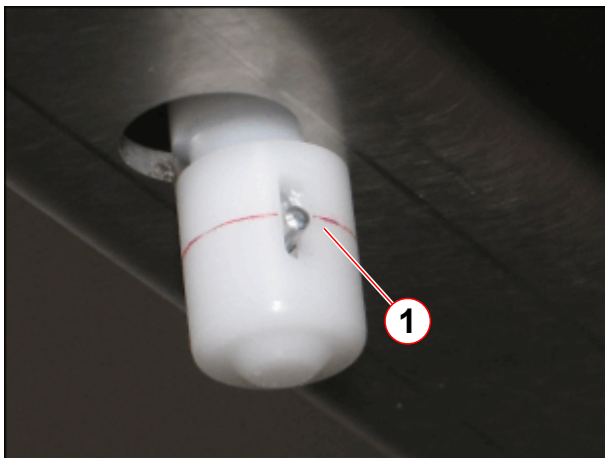

- (1) Correct tension of toothed belt (pin aligned with mark)

- Connect service aid to measure the tension of the toothed belt.
- Tighten the Allen screws (→ 1/Fig. 258 Page 156) until the pin is aligned with the mark (→ 1/Fig. 259 Page 156).
- When the toothed belt is adjusted to the proper tension, tighten the Allen screws (→ 2/Fig. 258 Page 156).
- Remove the service aid and place it on the store position (→ 1/Fig. 256 Page 155).

3.9.8 Final checks

- Switch ON F22 (RF-room) ptab and check for smooth horizontal movement of the tabletop plate.

The system should be on and SYNGO should be running.

- Move the table all the way out.
 - » The table should move constant and without any extraordinary noise.
 - » The table stops in the horizontal end position.
- Lower the table all the way.
 - » The table stops in the vertical end position.
 - » The table should move without any extraordinary noise.
- Move the table all the way up again.
 - » The table should move without any extraordinary noise.
 - » The table stops in the vertical end position.
- Move the table all the way in.
 - » The table should move constant and without any extraordinary noise.
 - » The table stops in the horizontal end position.
- Move the table to the Home position.

3.10 Toothed belt of horizontal drive

3.10.1 General

The "Partial drive" allows the horizontal table movement without whole body application in a limited movement range. A toothed belt is responsible for transferring the engine power into a horizontal table movement.

3.10.2 Tools and time required

3.10.2.1 Time

Approx. 60 min

3.10.2.2 Tools

Non-magnetic tool kit

Philips screwdriver

3.10.3 Preliminary

- Switch OFF F22 (RF-room) in the ACC.



According to the system/room specifications the movement range of the horizontal movement can be reduced.

- If the movement range of the horizontal movement is reduced it may be necessary to remove the end stops. (→ Fig. 289 Page 173)

3.10.4 Disassembly

3.10.4.1 Tabletop plate



Be careful and avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.



The tabletop plate may fall - support the tabletop plate

Basic table

Fig. 260: Coil plug cover (foot end)



- Remove the Allen screws on top of the cover from the coil plugs.


Fig. 261: Remove Allen screws



- Remove the lower Allen screws of the cover from the coil plugs.

Removable table


Fig. 262: RF-Cover RTT

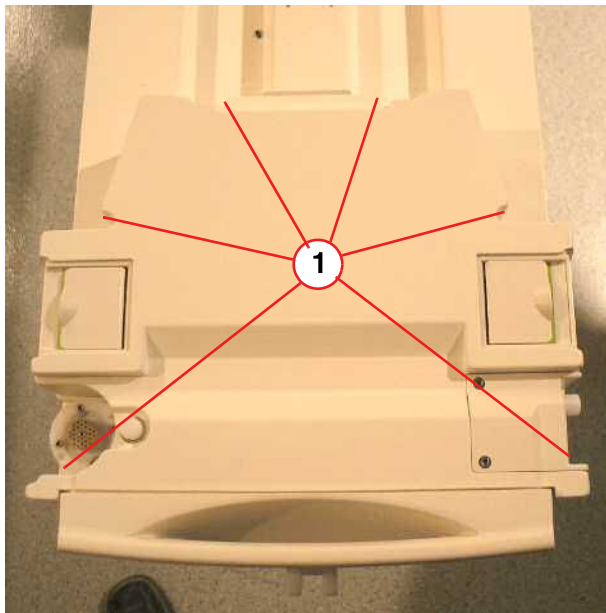

(1) Screws to be removed

- Remove the carrier plate of the RTT.
- Remove six screws on top of the cover from the coil plugs.
- Disconnect cables.

Fig. 263: Additional mounting bracket for foot end cover (newer systems)

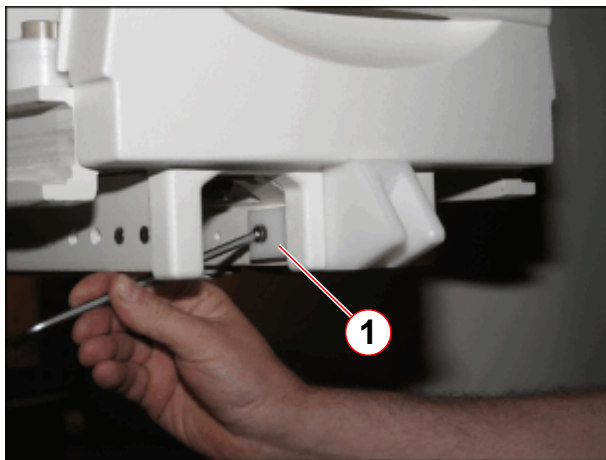

(1) Allen screw


- On newer systems, remove Allen screw of mounting bracket.

Fig. 264: Foot-end cover



- Remove screws of the foot-end cover.


Fig. 265: Foot end cover, removed

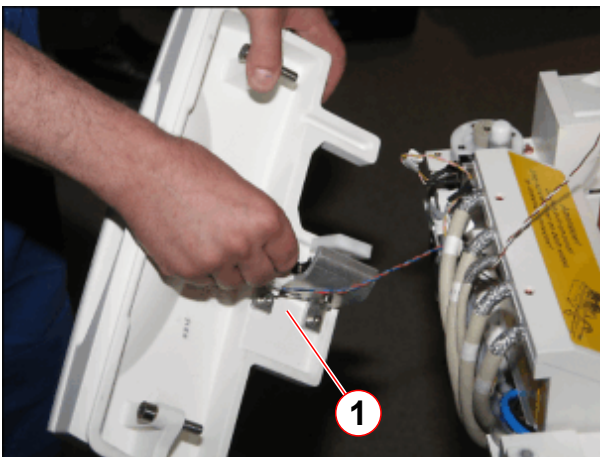

- Disconnect cable connections and remove the foot-end cover.
- Slide out the lower cover on the bottom of the patient table.

(1) Disconnect cables

Fig. 266: Mechanical end stop of partial drive unit

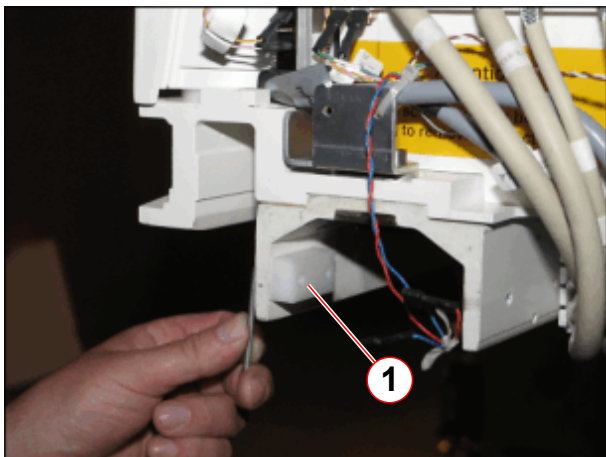

(1) End stop to be removed


- Remove mechanical end stop.

Fig. 267: Energy chain



- Remove the lower part of the energy chain.


Fig. 268: Tension unit for toothed belt

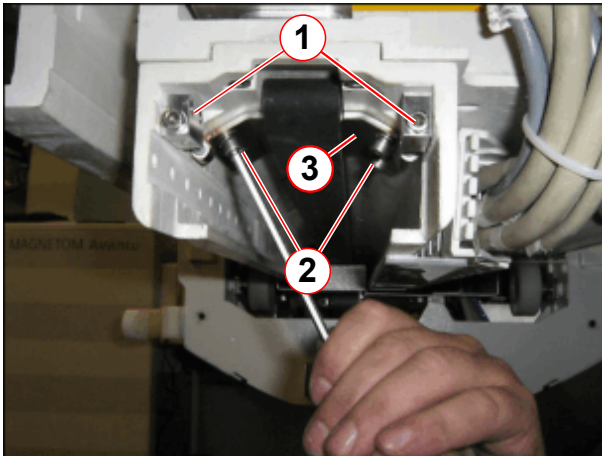

- (1) Allen screws for adjusting belt tension
- (2) Attachment screws
- (3) clamping plate

- Release tension of toothed belt, remove Allen screws and tooth belt tension unit.


Fig. 269: Cover of toothed belt mount (head end)

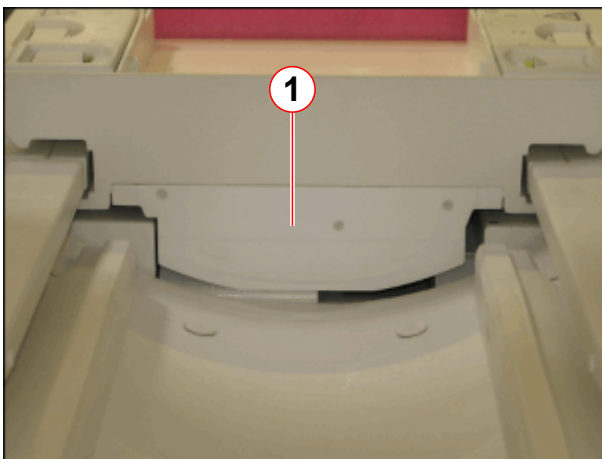

- (1) Cover of toothed belt mount

- Remove the cover of the toothed belt mount of toothed belt (1) on the head end from the ptab.


Fig. 270: Fixation of tooth belt - service end

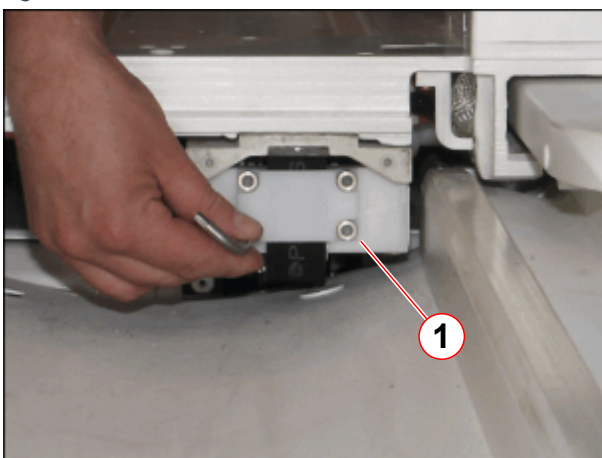

- (1) Allen screws of toothed belt mount

- Remove the mount of toothed belt on the head end from the ptab.

Fig. 271: Tabletop moved into the magnet

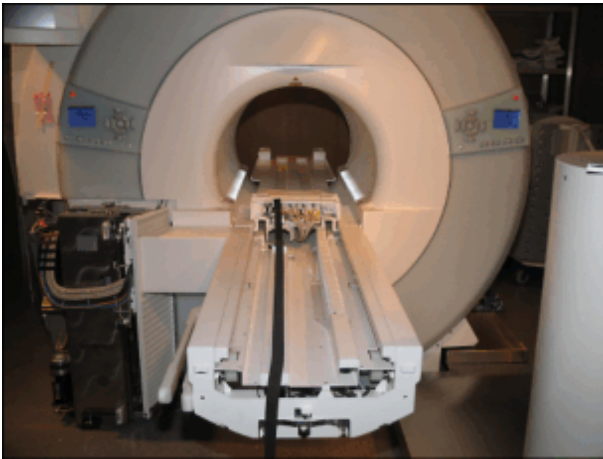


- Carefully push the tabletop plate manually into the magnet.

3.10.4.2 Covers / cables / energy chain



Be careful and avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.


Fig. 272: Plastic inner cover

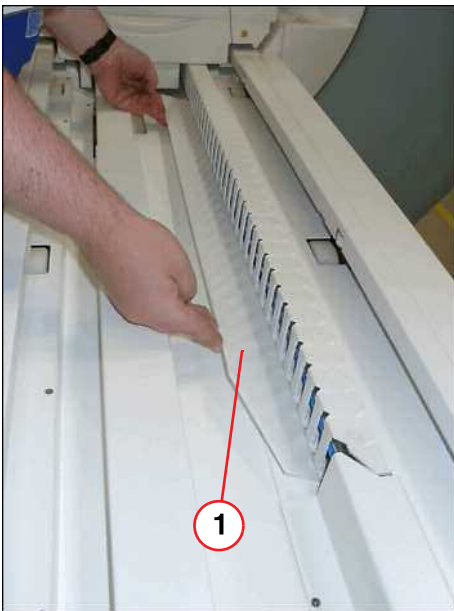

(1) Plastic inner cover

- Remove inner plastic cover.

Fig. 273: Remove cover of energy chain

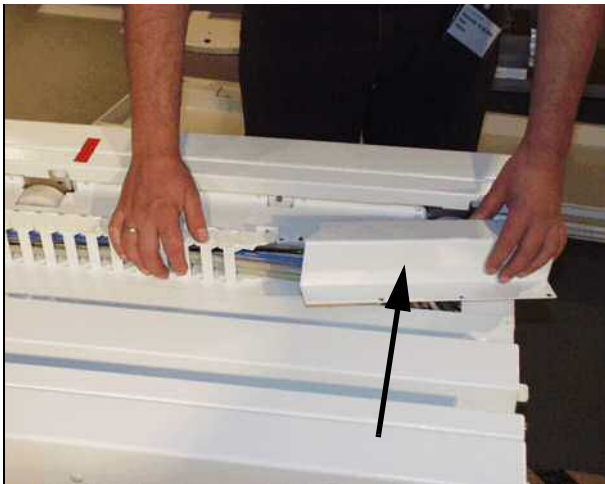


- Remove the upper energy chain cover.

Fig. 274: Remove energy chain




- Remove the upper energy chain of the bracket.

Fig. 275: Covers removed

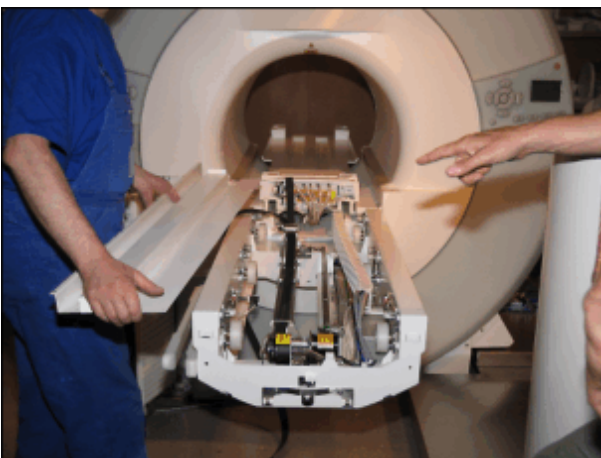


- Remove covers until the toothed belt is reachable.

Fig. 276: Toothed belt routing

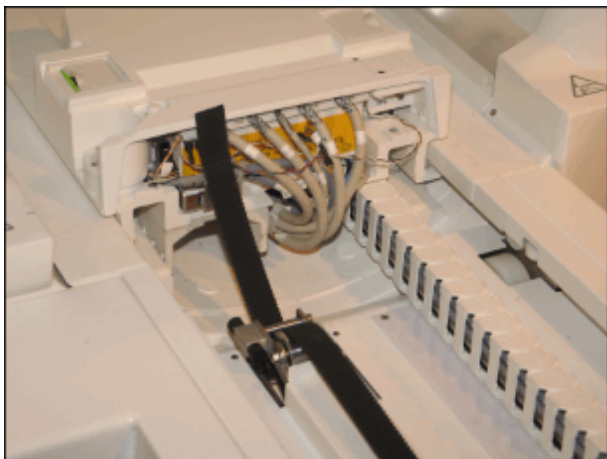


- Principle of toothed belt routing, head end table carrier.

Fig. 277: Schematic of the toothed belt routing

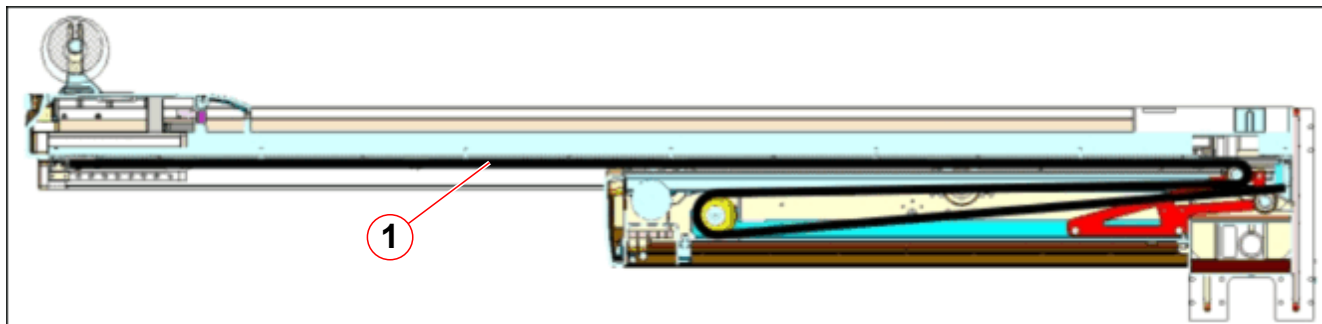

(1) Toothed belt

Fig. 278: Covers removed, ready for toothed belt exchange

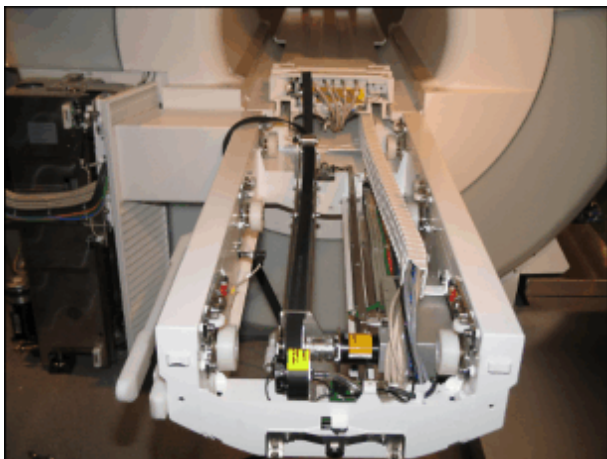


- Remove inner coves to get access to the toothed belt unit.
- Replace toothed belt.



Use the same length of screws when reinstalling the energy chain.


3.10.5 Assembly

Fig. 279: Tabletop into the magnet; toothed belt route

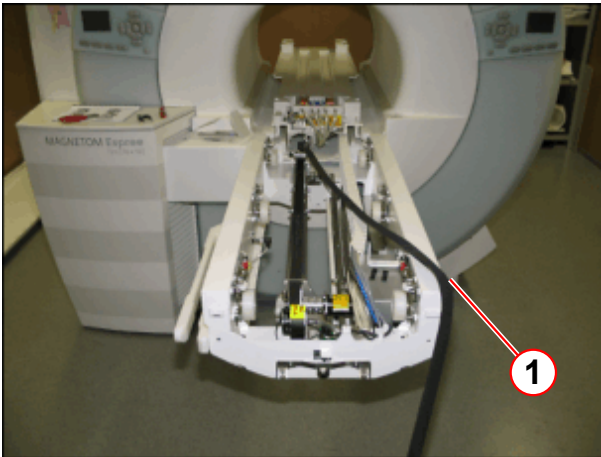

(1) To avoid misassembling, lay toothed belt over side


Fig. 280: Tabletop moved into the magnet

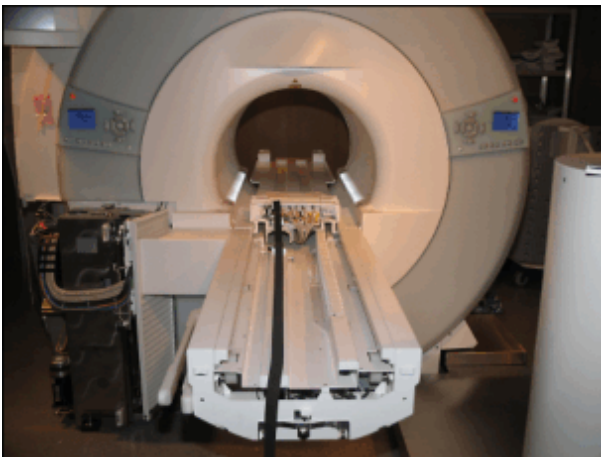

- Before you insert the the cover, ensure the correct toothed belt route.


- Make sure you remount the covers in the proper sequence and tighten all screws at the end of the assembly procedure.
- Install the energy chain.

Fig. 281: Toothed belt attachment (head end)

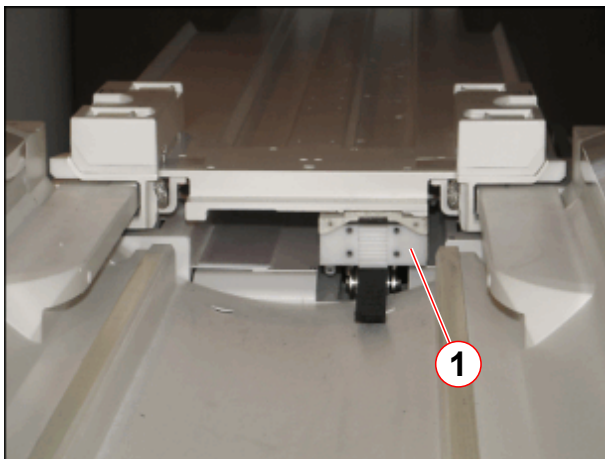


(1) Toothed belt mount (head end)

Fig. 282: Fixation of tooth belt - service end

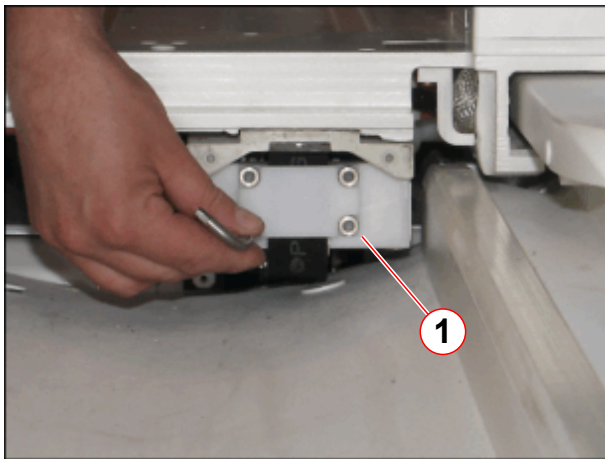

(1) Allen screws of toothed belt mount

Fig. 283: Cover of toothed belt mount (head end)

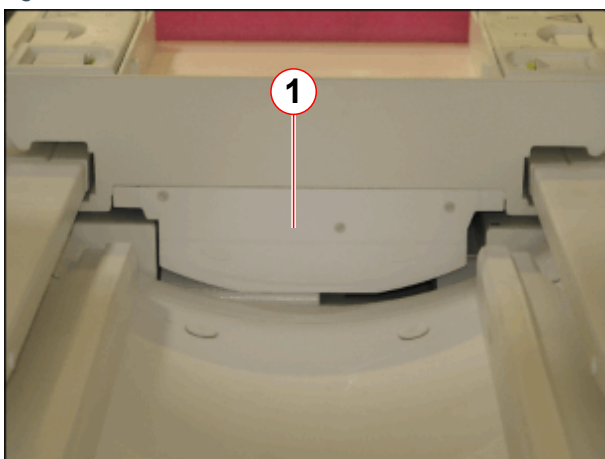

(1) Cover of toothed belt mount

- Carefully pull the tabletop plate manually back to the food end side (patient end).
- Remount the mount of toothed belt on the head end from the ptab.


- Install the cover of the toothed belt mount (1) on the head end from the ptab.

Fig. 284: Mounting the toothed belt (foot end)



- Principle of toothed belt tension unit, located on the foot end side of the ptab.

Fig. 285: Service tool for toothed belt tension

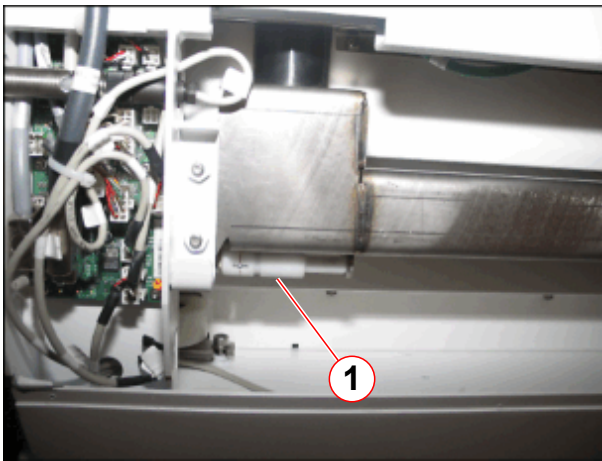

- The service aid for adjustment of the toothed belt tension is located underneath the table top.
- Snap out the service aid for adjustment of the toothed belt tension.


(1) Storage position of service tool for toothed belt tension

3.10.6 Adjustments

Fig. 286: Service tool, toothed belt tension




Fig. 287: Tension unit for toothed belt

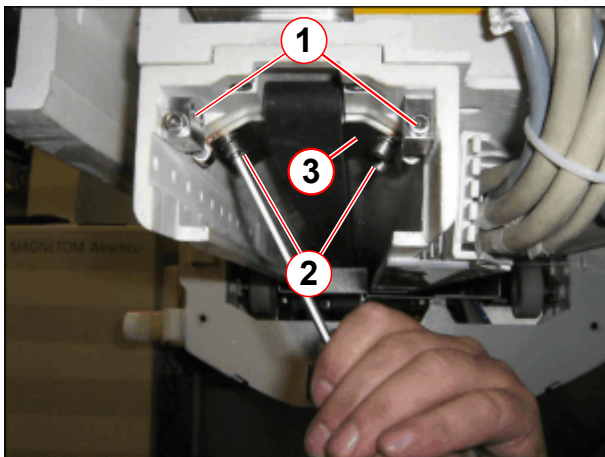

- (1) Allen screws for adjusting belt tension
- (2) Attachment screws
- (3) clamping plate

Fig. 288: Service tool for toothed belt tension in use

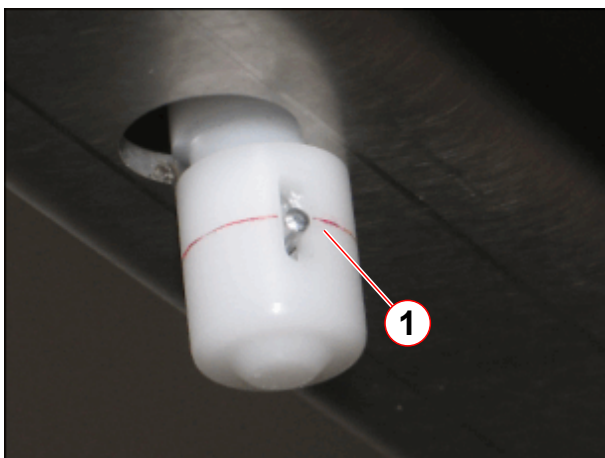

- (1) Correct tension of toothed belt (pin aligned with mark)

- Connect service aid to measure the tension of the toothed belt.
- Tighten the Allen screws (→ 1/Fig. 287 Page 170) until the pin is aligned with the mark. (→ 1/Fig. 288 Page 170)
- When the toothed belt is adjusted to the proper tension, tighten the Allen screws. (→ 2/Fig. 287 Page 170)
- Remove the service aid and place it on the store position. (→ 1/Fig. 285 Page 169)

3.10.7 Concluding steps

- Install the end stops if necessary.
- Install the remaining covers.
- Install foot end plate and the cover of the coil plugs.

3.10.8 Final checks

- Please refer to final checks horizontal drive unit. (→ Final checks / Page 157)

3.11 Sound transducer

3.11.1 General

The sound transducer is responsible for patient communication. The sound transducer can be replaced. Consider the cable run during assembling.



CAUTION

The sound transducer consists of magnetic objects!

Magnetic objects may cause injury.

» **Maintain a safe distance to the magnetic field!**

3.11.2 Tools and time required

3.11.2.1 Time

Approx. 60 min

3.11.2.2 Tools

Non-magnetic tool kit

Philips screwdriver, midsize

3.11.3 Preliminary

- Switch OFF F22 (RF-room) in the ACC.



According to the system/room specifications the movement range of the horizontal movement can be reduced.

With the whole body version of the horizontal drive, the movement range of the horizontal movement may be reduced. If that is the case, remove the end stops.

(→ Fig. 289 Page 173)

Before disassembly the patient table covers on a system with a partial drive unit: The toothed belt must be released to push the table top to the service end of the magnet. Therefore refer to (→ Prerequisite for cover removal / partial drive / Page 98).

3.11.4 Disassembly

3.11.4.1 Covers



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.


- For disassembling/assembling the patient table covers we recommend a mid-sized Philips screwdriver.
- Push the tabletop plate manually into the magnet
- If necessary, remove the mechanical end stops

Fig. 289: Table carrier underside (whole body, limited movement range)




(1) Mechanical end stops

Fig. 290: Foot-end cover

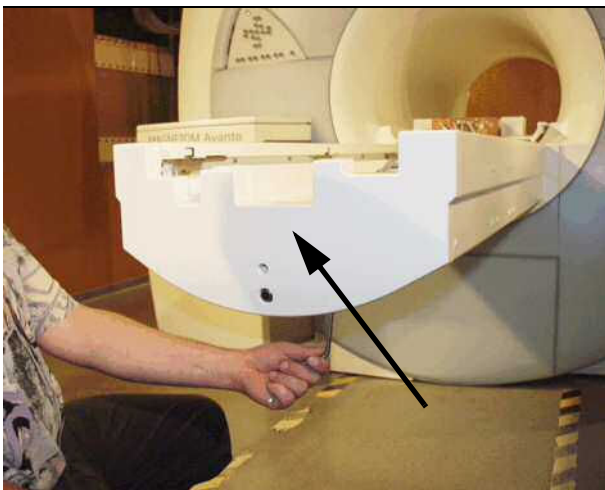

- Remove the foot-end cover of the patient table.


Fig. 291: Plastic inner cover

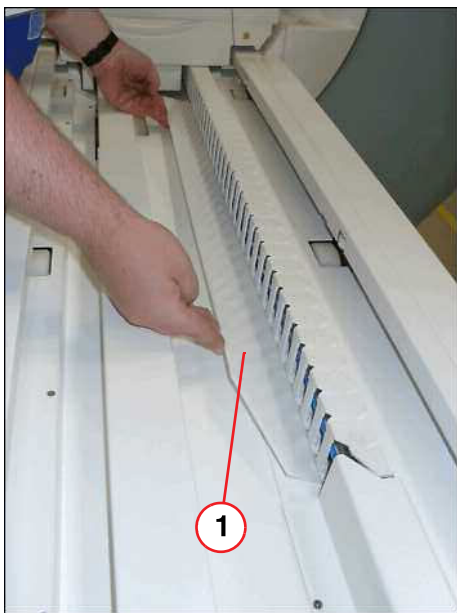

(1) Plastic inner cover

- Remove inner plastic cover.

Fig. 292: Remove cover of energy chain

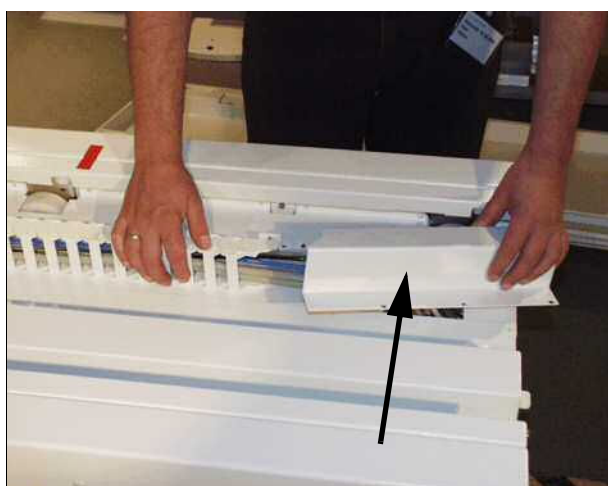


- Remove the upper energy chain cover.

Fig. 293: Remove energy chain



- Remove the upper energy chain of the bracket.


Fig. 294: Inner main cover

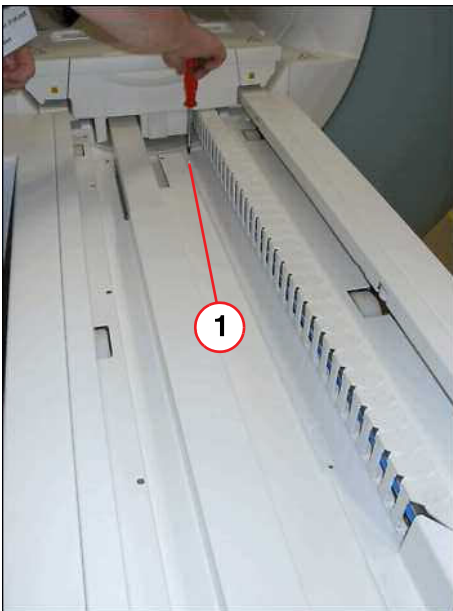

- Remove screws of main inner cover.

(1) Screws to be removed

Fig. 295: Allen screws

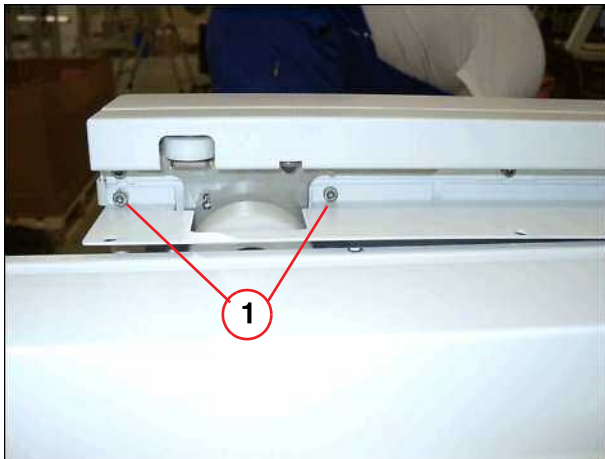


(1) Screws to be removed

- Remove two Allen screws at the left side of the table end (foot) before you remove the main inner cover.


- Remove main inner cover.

Fig. 296: Inner main cover

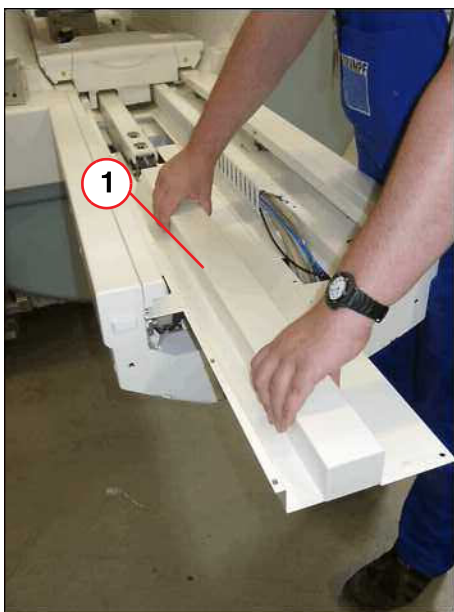

(1) Inner main cover

Fig. 297: Button cover

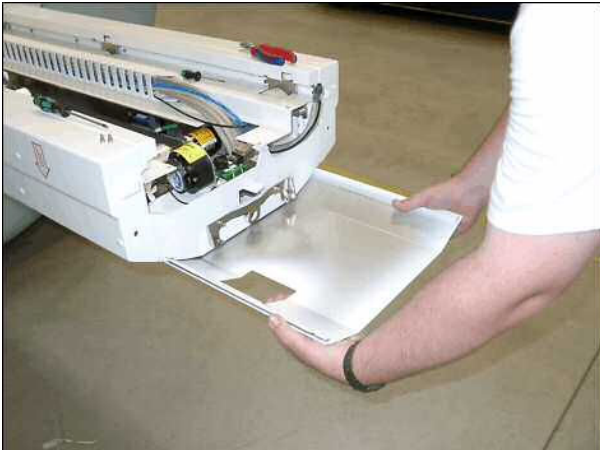

- Remove bottom cover.

3.11.4.2 Sound transducer


Fig. 298: Cut air hose



- Cut black air hose approx. 5 cm.

- Remove four screws, air hose and sound transducer.

Fig. 299: Sound transducer

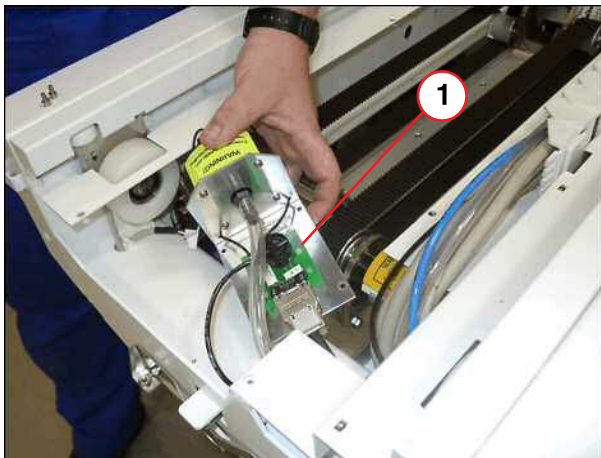

(1) Sound transducer to be removed



CAUTION


The sound transducer consists of magnetic objects!
Magnetic objects may cause injury.

» Maintain a safe distance to the magnetic field!

3.11.5 Assembly

3.11.5.1 Sound transducer


Fig. 300: Sound transducer

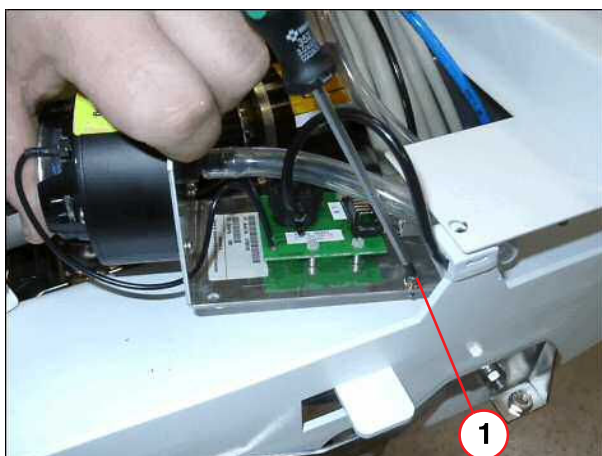

(1) Insert first Allen screw

- Insert first one Allen screw as shown.

Fig. 301: Cable run

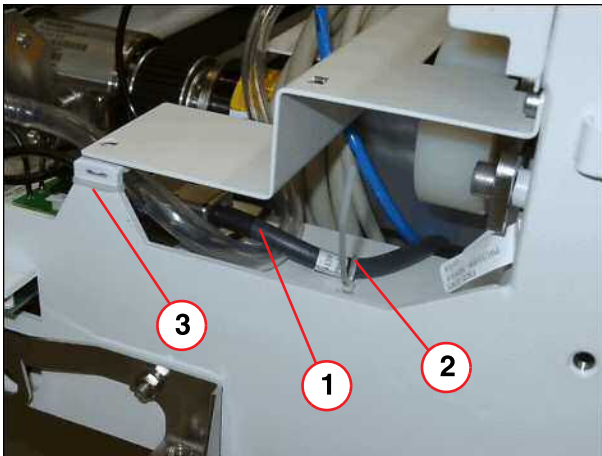

- (1) Run the black cable above the air hose
- (2) Secure black cable with tie wrap
- (3) Fix small black air hose with tie wrap

- To prevent the cable from being crushed, pay attention and route the cables as shown, e.g. the black cable x 23 above the big air hose. Safety small black air hose and black cable x 23 as shown.

3.11.6 Concluding steps

For assembling the partial drive, please refer to (→ Assembly / Page 167), not valid for Trio A Tim Systems.

- Reinstall the covers in reverse order, tighten all screws at the end of the assembling procedure.
- If necessary, reinstall the mechanical end stops.
- Switch ON F22 (RF-room).

3.11.7 Adjustments

For adjustment of the toothed belt, please refer to (→ Adjustments / Page 170), not valid for Trio A Tim Systems.

3.11.8 Final checks

- Check for smooth patient table movement.
- Do a check of the function of the patient communication.

3.12 RF coil plugs

3.12.1 General

The coil plugs at the foot (magnet front) and on the head (magnet end) positions of the patient table can be replaced.

The coil plugs are covered with caps, each of which can be replaced in case of damage.

3.12.2 Tools and time required

3.12.2.1 Time

Approx. 30 min/ 30 min

3.12.2.2 Tools

Non-magnetic tool kit

Service plug

3.12.3 Preliminary

- Switch OFF F24 (RFIS)
- Switch OFF F22 (RF-room) in the ACC

3.12.4 Disassembly

3.12.4.1 Coil plugs at the foot end side (magnet front)



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

- If the movement range of the horizontal movement is reduced it may be necessary to remove the end stops. (→ Fig. 289 Page 173)

Basic table

Fig. 302: Coil plug cover (foot end)



- Remove the Allen screws on top of the cover from the coil plugs.


Fig. 303: Remove Allen screws



- Remove the lower Allen screws of the cover from the coil plugs.

Removable table

Fig. 304: RF-Cover RTT

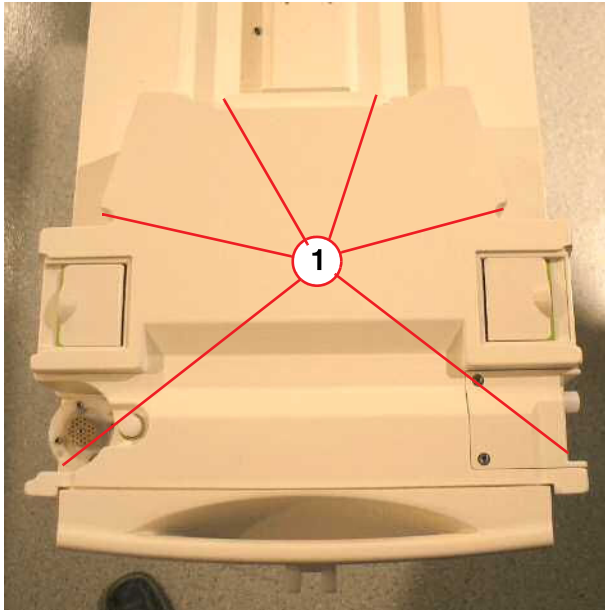

(1) Screws to be removed

- Remove the carrier plate of the RTT.
- Remove six screws on top of the cover from the coil plugs.
- Disconnect cables.

Fig. 305: Foot-end cover



- Remove the two M8 Allen screws, if necessary disconnect cables, and remove foot end plate.

3.12.4.2 RF coil plugs



Observe the information in the note when disassembly the coil plugs.



Do not pull on cord, use a flat screwdriver on both sides to remove connector.

Fig. 306: RF cover from interconnect box



- Remove the RF cover from the interconnect box.

- Loosen the screws of defective RF coil plugs.
- Coil plug 7-8 belongs to X101.
- Coil plug 9-10 belongs to X104.
- Release the corresponding strain relief from the coil plug (X101 or X104).

- Disconnect RF cables, use a flat screwdriver to open the RF connectors.

Fig. 307: Cable connections



- Install the corresponding covers.

3.12.5 Assembly



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

- Replace the defective RF coil plugs including the RF wave trap.
- Install the new RF coil plug assembly.
- Reconnect RF cables.



Be careful to avoid pinching any cables when assembling the cover of the coil plug connectors.

- Make sure the RF cable connections are installed properly.
- Reinstall the strain relief from the replaced coil plug.

3.12.6 Concluding steps

- Reinstall the RF cover of the interconnect box.

- Install the corresponding covers.

3.12.7 Final checks

- Switch ON F24 (RFIS).
- Start QA coil check and do a check of the replaced coil plug.

3.12.8 Disassembly

3.12.8.1 Coil plugs at the head end (back of the magnet) / Coil plug caps



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.



Do not pull on cord, use a flat screwdriver on both sides to remove the connector

- Manually push the tabletop plate into the magnet until the head end of the tabletop plate is accessible.

- Remove the screws of the affected coil plug.

Fig. 308: Screws of coil plugs (magnet end, back view)



- Remove the cable strain relief of the affected coil plug.

Fig. 309: Strain relief of coil plug cable (magnet end, back view)

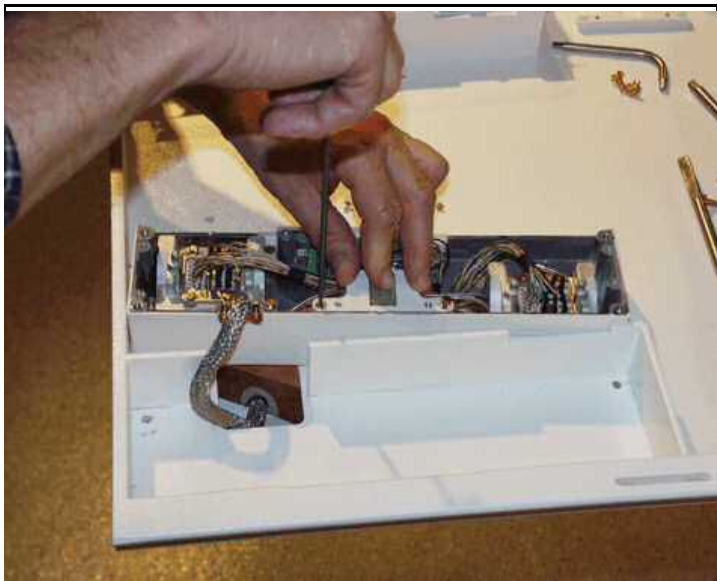


Fig. 310: Coil plug assembly (head end)



- Remove the affected coil plug assembly.

Fig. 311: Plug fastener of coil plugs (magnet end)




- Remove the RF cover of the affected coil plug assembly by loosening the screws.

3.12.8.2 Replacing the coil plug caps

The coil plugs are covered with caps, each of which can be replaced in case of damage.

- Remove inner cover of defective coil plug.


Fig. 312: Remove cover

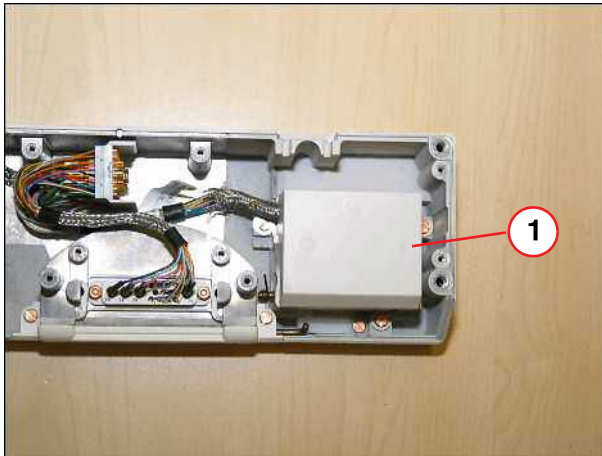

(1) Cover of coil plug assembly


- Remove the screws of the defective coil plug cap.

Fig. 313: Coil plug housing

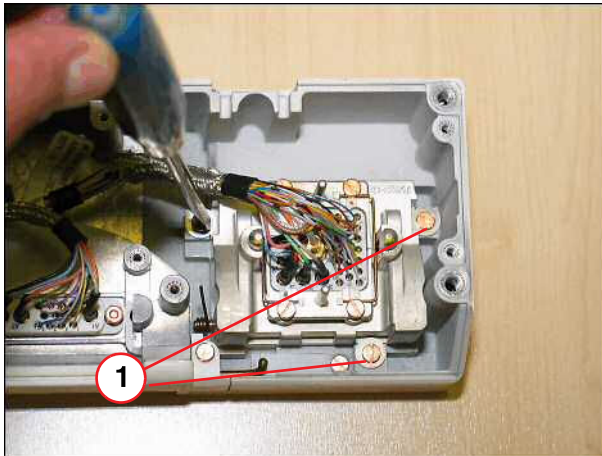

(1) Screws to be removed


Fig. 314: Coil plug

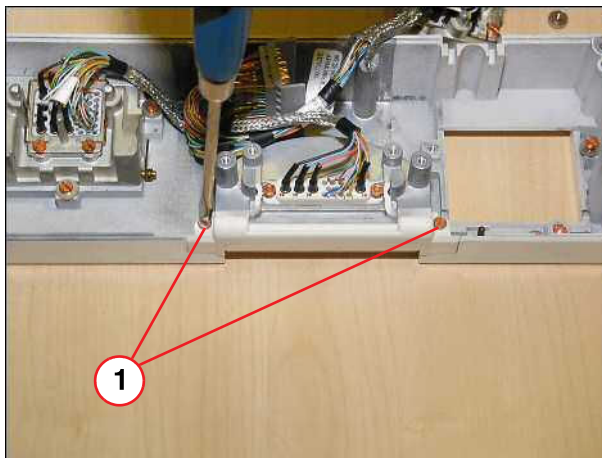

(1) Screws to be removed

- Replace the defective coil plug cap.

Fig. 315: Coil plug

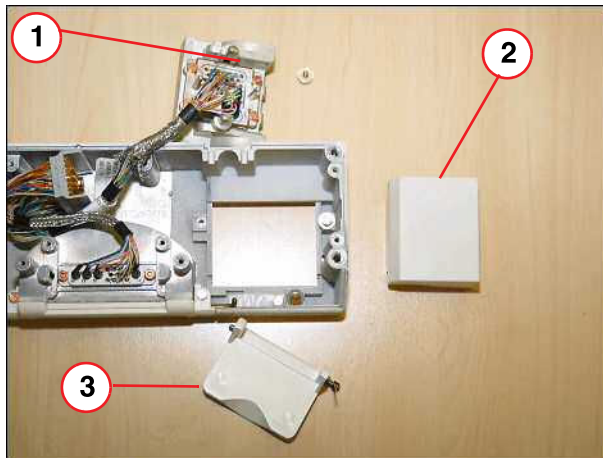

- (1) Coil plug housing
- (2) Cover of coil plug housing
- (3) Coil plug cap

- Remove the coil plug fastener.
- Disconnect RF cables; use a flat screwdriver to open the RF connectors.

Fig. 316: Cable connections (magnet end)



3.12.9 Assembly

- Replace the defective coil plug cap.
- Replace the defective coil plug assembly.
- Replace protective film (part of the coil plug assembly delivery) if necessary. For additional information, refer to UI 26/04.
- Connect the RF cables and install the RF connector.
- Install the coil plug fastener.



CAUTION

Bend the cable as shown in the picture below.

Any contact between the shield of the cable and the housing of the RF trap must be avoided.

» Otherwise an overheating could happen in the components.

Fig. 317: RF traps

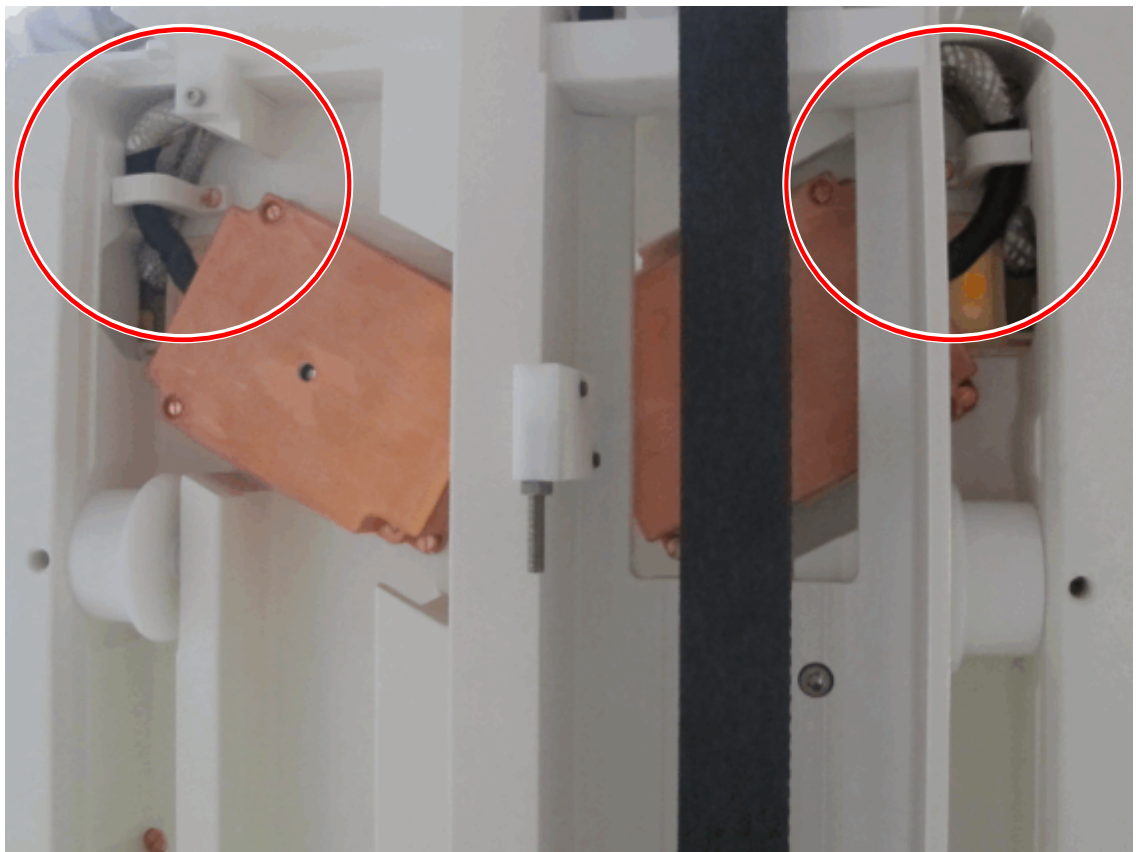


Fig. 318: Correct cable routing, left RF trap



Fig. 319: Correct cable routing, right RF trap



3.12.10 Concluding steps

- Reinstall the RF cover.



Be careful to avoid pinching any cables when assembling the cover of the coil plug connectors.

- Reassemble the coil plug assembly.
- Install the coil plug assembly to the tabletop plate.

3.12.11 Final checks

- Switch ON F24, F22.
- Start QA coil check .

3.13 Control switch S30

3.13.1 General

The S30 maintains the vertical collision identification. If the S30 detects uncoupling of the spindle nut, the switching contact is open. In case the S30 has to be replaced, replacement is described in this part.

3.13.2 Tools and time required

3.13.2.1 Time

Approx. 60 min

3.13.2.2 Tools

Non-magnetic tool kit

3.13.3 Preliminary

- Switch OFF F22 in the ACC.

3.13.3.1 Location and principle of S30

- The switch location of the S30 is shown in the following figures.

Fig. 320: Princip of upper spindle part


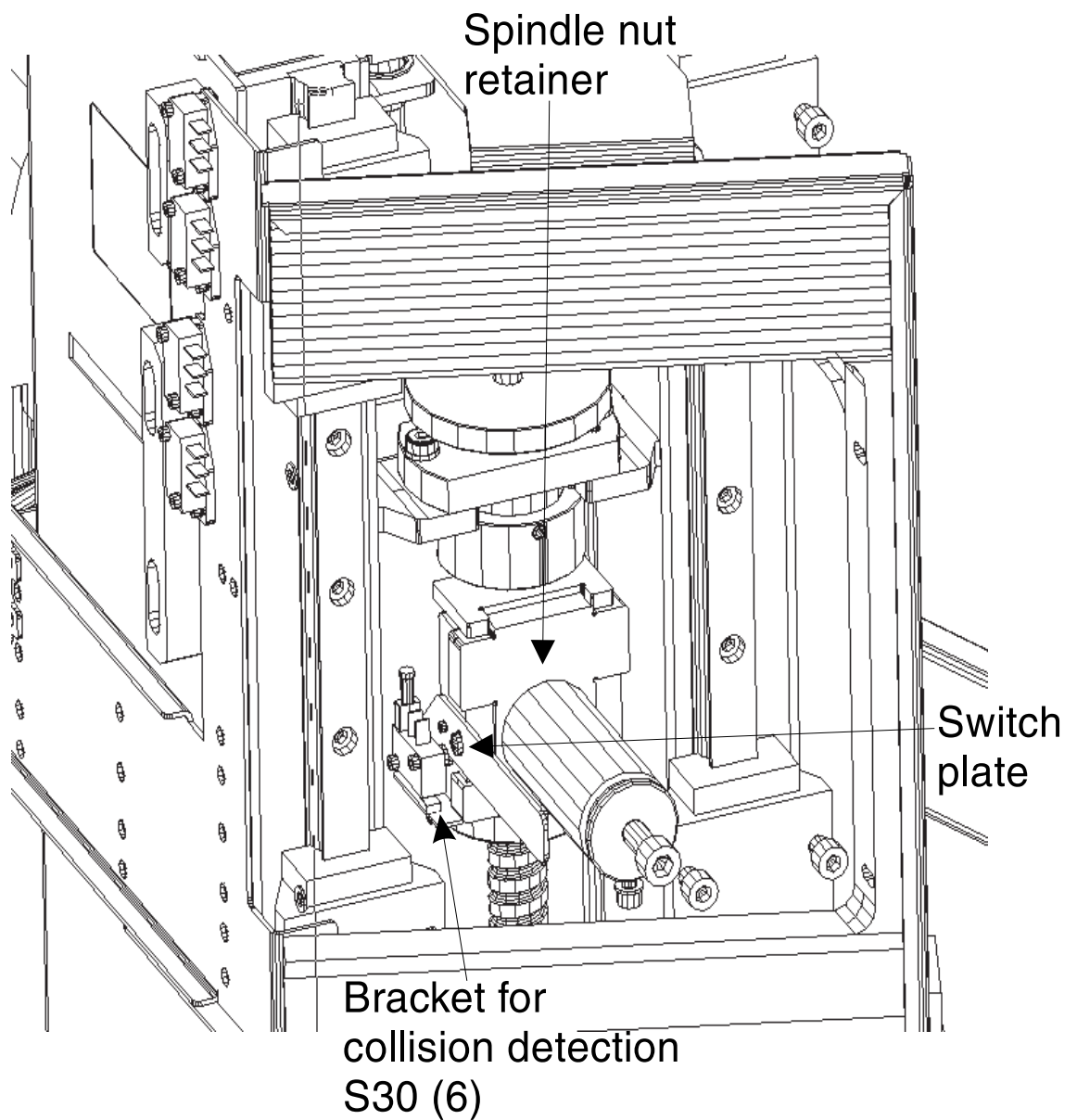



Fig. 321: Spindle with S30

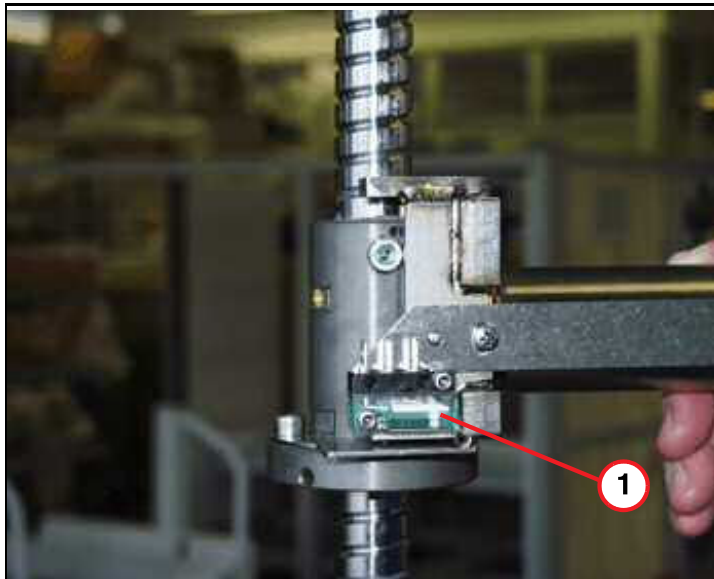

(1) Principle of S30, normally closed

- The principle of S30, normal state= switch is actuated, switching contact is closed, rotate mounting bracket for switch adjustment.

Fig. 322: S30 assembly

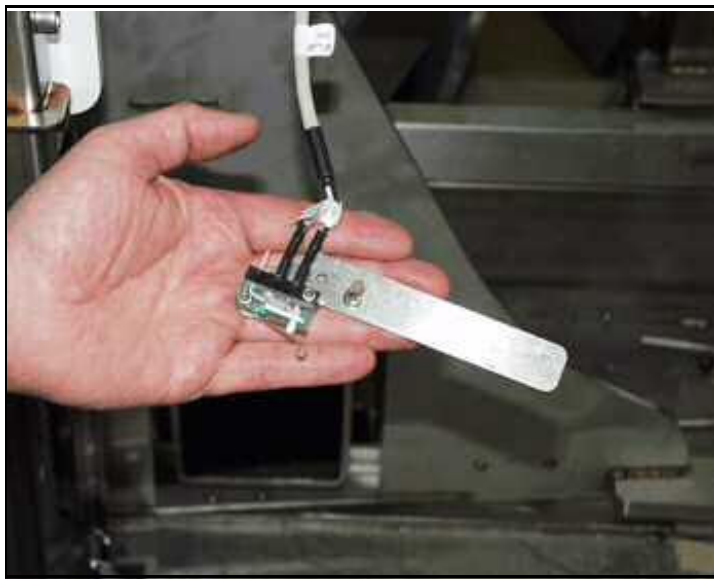

- S30 assembly.

3.13.4 Disassembly

3.13.4.1 Covers



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

The removal of the service access cover and the lifting column cover is described in the REP/system, refer to (→ System / M6-020.841.10).

Fig. 323: Lower magnet, front cover



- Remove lower front cover of the magnet.

3.13.4.2 Bellows guide

Fig. 324: Remove the M6 Allen screws of the bellows guide



- Separate the bellows from the adhesive tape and remove the M6 Allen screw from the bellows guide.
- Remove the M6 Allen screws of the bellows guide.

Fig. 325: Lower Allen screws of bellows guide



- Remove the lower M5 Allen screws.

Fig. 326: Remove hex bolt of bellows guide




- Remove the M6 Hex bolt of the bellows guide
- Separate the bellows from the adhesive tape and remove the M6 Allen screw from the bellows guide.
- Remove the bellows guide after removing all screws.

3.13.4.3 Switch assembly S30

- If necessary cut off tie wraps from the cables of the lifting column.

- Beginning with patient table serial number 1178, the lifting column is equipped with a service cut-out to access the S30 switch.


Fig. 327: Cut-out on lifting column

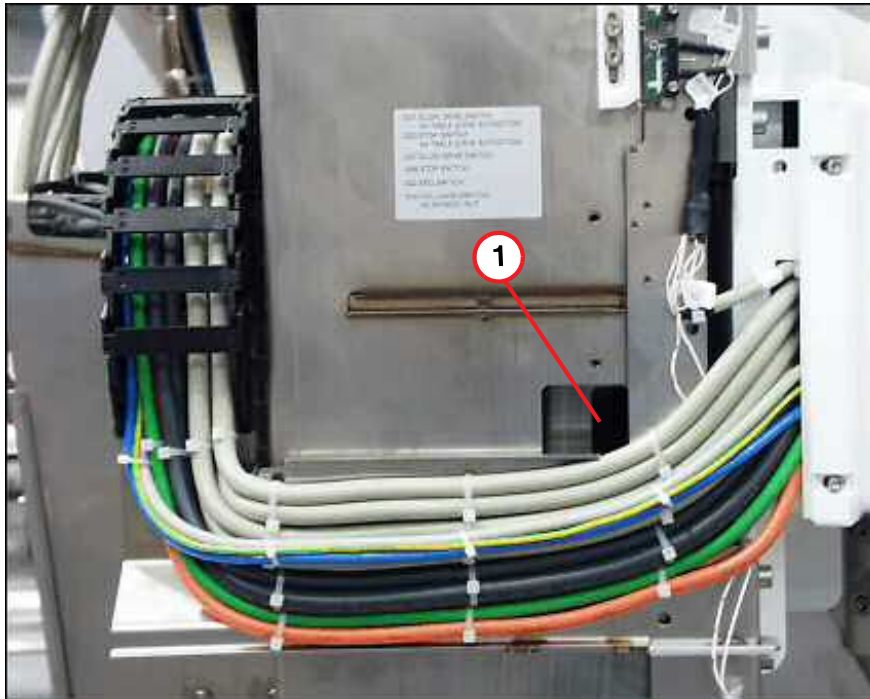

(1) Access to S30


- Use a Phillips screwdriver to remove the S30 assembly.

Fig. 328: Lifting column, location of S30

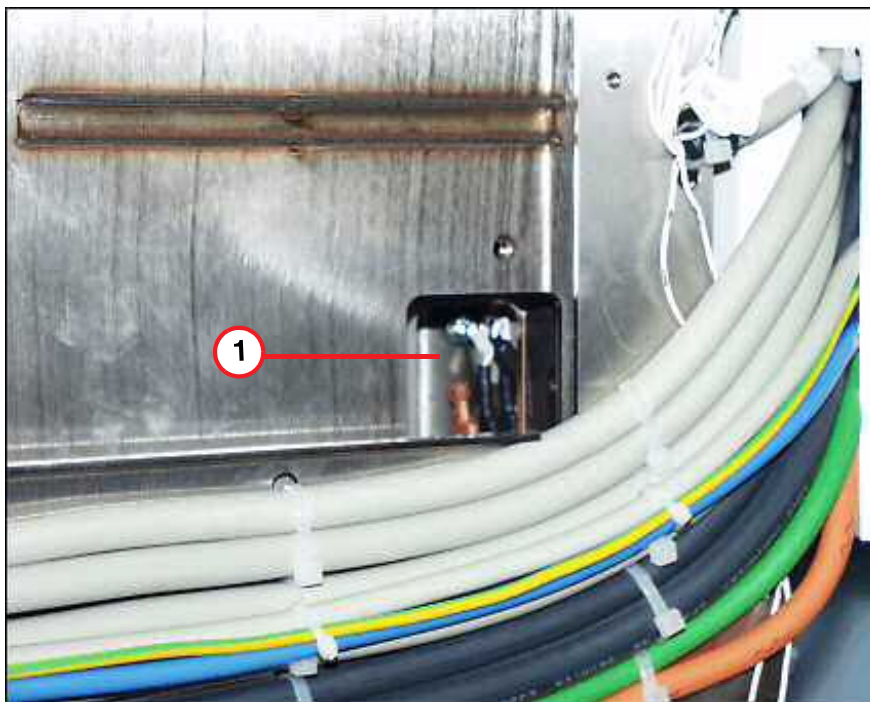

(1) S30


Fig. 329: S30 attachment

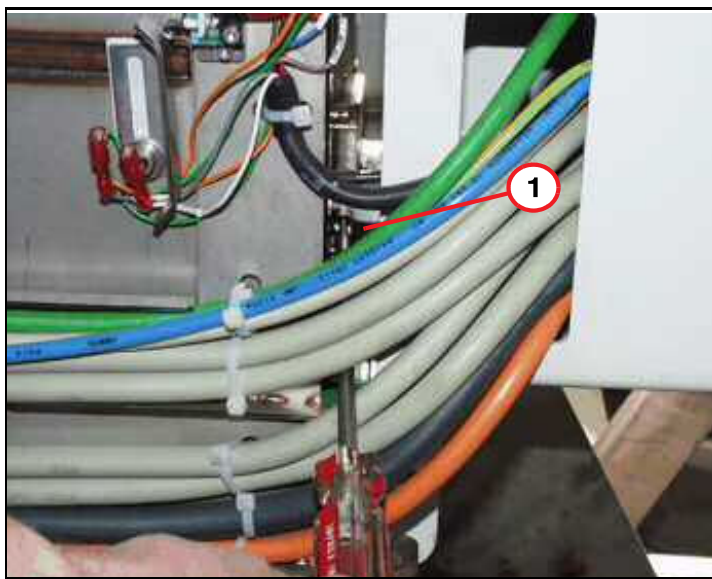

(1) Feed the Phillips screwdriver through the cables

- For older systems use a Philips screwdriver to remove the S30 assembly.

Fig. 330: S30 mounting bracket

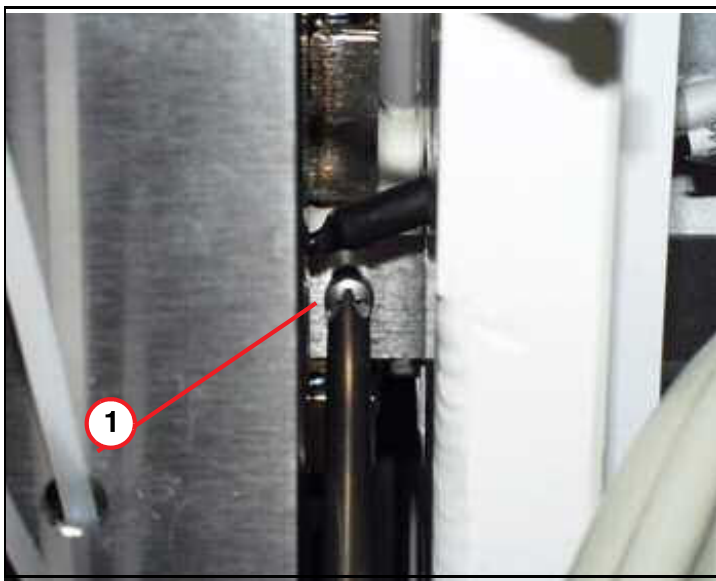

(1) Tighten screw after adjustment

Fig. 331: S30 switch assembly



- Catch the S30 switch assembly from underneath the support beam and pull it out.

3.13.5 Assembly




Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

Fig. 332: S30 switch assembly



- Disconnect cables and replace the defective S30 switch assembly.
- Feed in the new S30 switch assembly.


Fig. 333: Location of S30

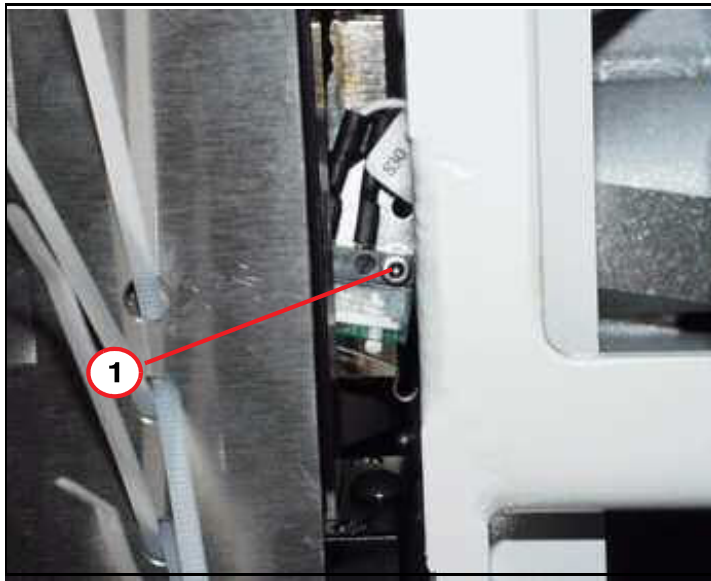

(1) Mounting screw and rotation point of S30

- Slide the S30 switch assembly through the support beam to the spindle nut for fastening.

Fig. 334: S30 mounting bracket

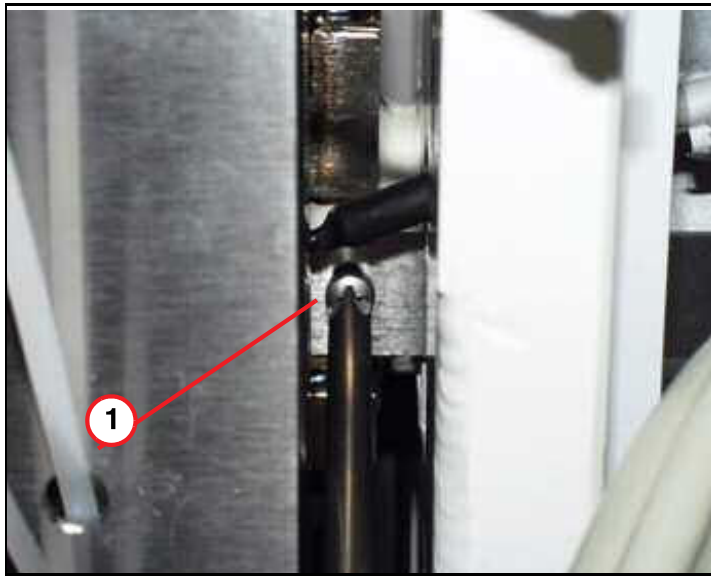

(1) Tighten screw after adjustment

- Safety the S30 switch assembly with the Philips screw and adjust the S30 switch.
- If the switch assembly is properly adjusted, S30 is always closed.

3.13.6 Concluding steps

- If necessary reconnect cables and install cables to the lifting column.
- Reinstall the covers.

3.13.7 Final checks

- Switch ON F24 (RFIS).

- Start patient table movement and check for smooth operation.
- Start patient table service menu for switch test and check S30 HI status.
- Start to move the table down until it reaches the end stop.
- Insert hand crank and lower the table some more until the switch status of the S30 changes from HI to LOW in the Service mode.

3.14 Cable replacement

3.14.1 General

The replacement of a defective cable, running in the energy chain through the patient table, is described in this part.

3.14.2 Tools and time required

3.14.2.1 Time

Approx. 180 min

3.14.2.2 Tools

Non-magnetic tool kit

Tie wraps

3.14.3 Preliminary

- Switch OFF the corresponding circuit breaker.
- Switch OFF F22 (RF-room) in the ACC.



According to the system/room specifications the movement range of the horizontal movement can be reduced.

With the whole body version of the horizontal drive, the movement range of the horizontal movement may be reduced. If that is the case, remove the end stops.
(→ Fig. 289 Page 173)

Before disassembly the patient table covers on a system with a partial drive unit: The toothed belt must be released to push the table top to the service end of the magnet. Therefore refer to (→ Prerequisite for cover removal / partial drive / Page 98).

3.14.4 Disassembly

3.14.4.1 Covers



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

- Remove the service access cover and the lifting column cover


The removal of the service access cover and the lifting column cover is described in the REP/system, refer to (→ System / M6-020.841.10).

Fig. 335: Lower magnet, front cover



- Remove lower front cover of the magnet.

Fig. 336: Foot-end cover

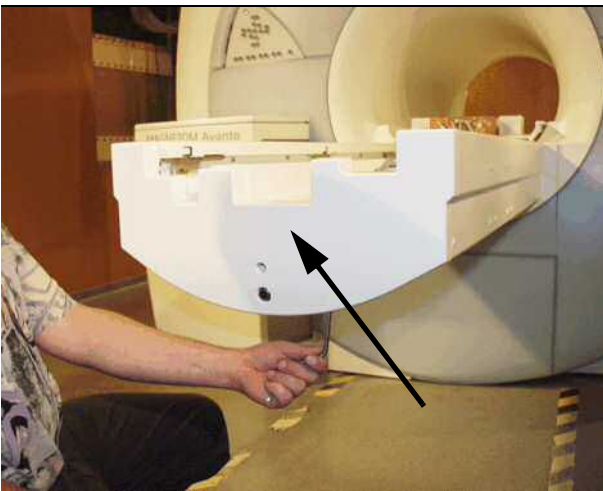

- Remove the foot-end cover of the patient table.

Fig. 337: Remove Allen screws



- Remove the lower Allen screws of the cover from the coil plugs.

Fig. 338: Coil plug cover (foot end)



- Remove the Allen screws on top from the cover of the coil plugs.

3.14.4.2 Bellows guide

Fig. 339: Remove the M6 Allen screws of the bellows guide



- Separate the bellows from the adhesive tape and remove the M6 Allen screw from the bellows guide.
- Remove the M6 Allen screws of the bellows guide.

Fig. 340: Lower Allen screws of bellows guide



- Remove the lower M5 Allen screws.

Fig. 341: Remove hex bolt of bellows guide




- Remove the M6 Hex bolt of the bellows guide.
- Separate the bellows from the adhesive tape and remove the M6 Allen screw from the bellows guide.
- Remove the bellows guide after removing all screws.

3.14.4.3 Cables/energy chain



Be careful to avoid damaging cables while removing the covers of the patient table. Use caution when moving the patient table.

Fig. 342: Remove cover of energy chain

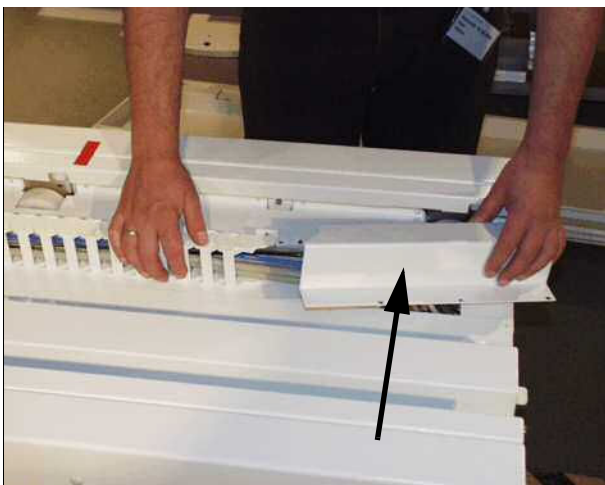


- Remove the upper energy chain cover.

Fig. 343: Remove energy chain



- Remove the upper energy chain of the bracket.

Fig. 344: Inner cover

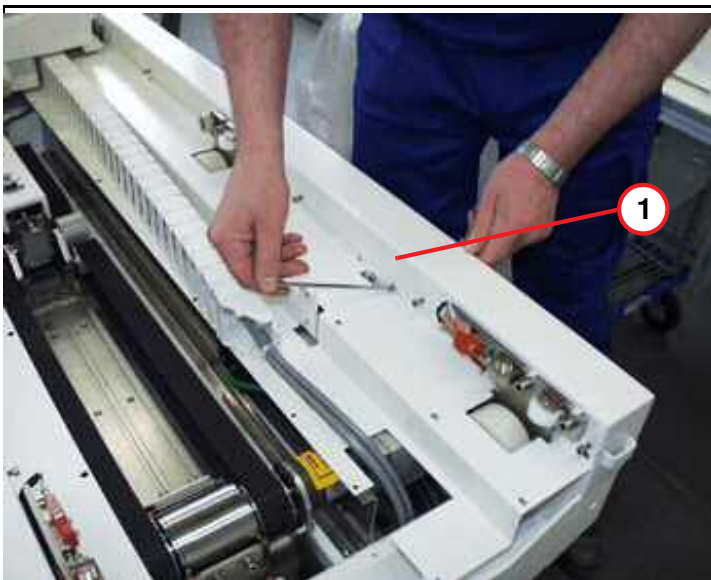

- Remove the right guide rail cover and the inner top cover.

(1) Top inner cover


- Remove the inner main cover.

Fig. 345: Inner cover



Fig. 346: Inner cover

- Remove the inner right cover.

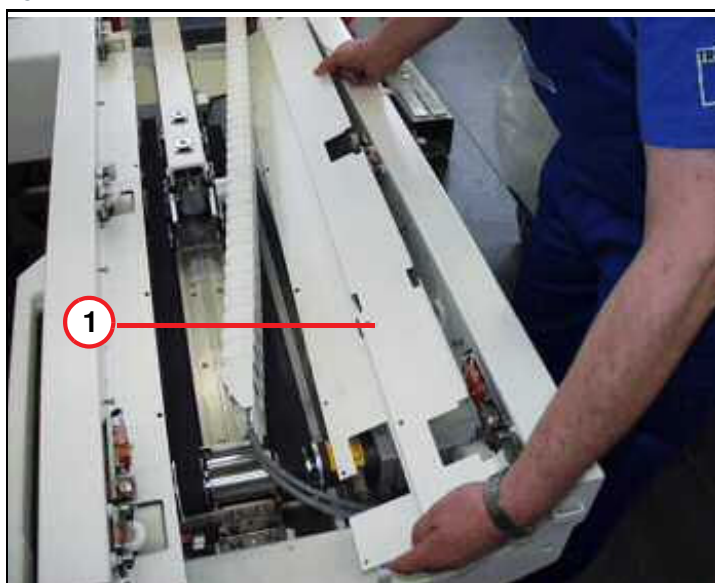


(1) Inner cover, right side

Fig. 347: Energy chain for lifting column



- Open chain links from the energy chain on the lifting column. If necessary cut tie wraps and remove the energy chain. Also, remove the cover of the patient table electronics.

Fig. 348: End piece of the energy chain

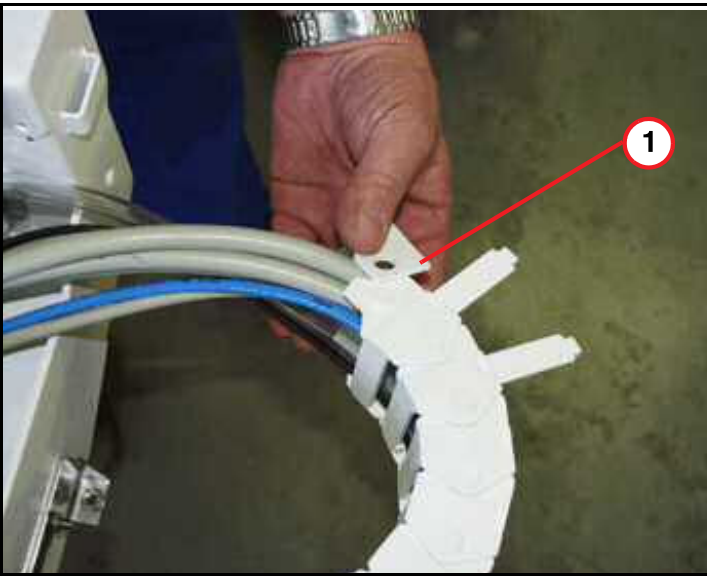


- Disconnect defective cable from the corresponding connector and release defective cable on both sides.
- Remove end piece of energy chain from both sides at the table.

(1) Remove end piece from both sides



The replaced cable has to have the same length as the old one. Therefore measure before replacing it. Measure from end to end, that is, from the end of the energy chain to the end of the affected cable at table side.

Fig. 349: Chain links of the energy chain

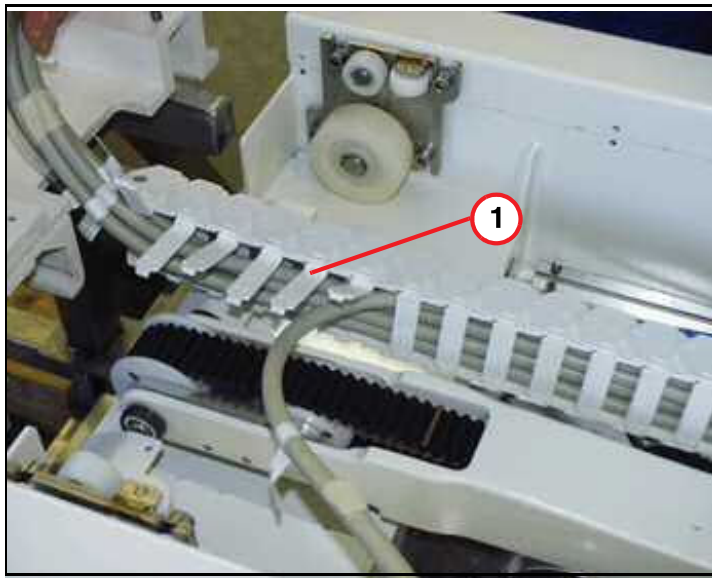

(1) Open chain links of the energy chain

- Open the energy chain and remove the defective cable from the chain at table side.

3.14.4.4 Assembly

For assembling the partial drive, please refer to (→ Assembly / Page 167), not valid for Trio A Tim Systems.




The replaced cable must be the same length as the old one. Therefore, measure it before replacing it. Measure from end to end; i.e., from the end of the energy chain to the end of the affected cable at the side of the table.



To avoid damage on the cables in the moving parts (energy chain), the replaced cable should run parallel with the other cables.


Fig. 350: Cable replacement

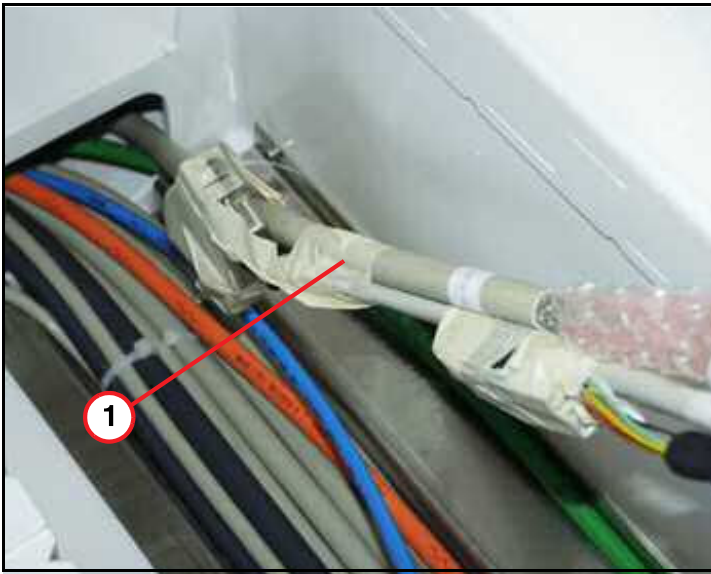

(1) Tape the new cable to the defective one

- Tape the new cable to the defective one at table side. If necessary cut off connector of defective cable.
- From the magnet side, pull the new cable through the tabletop plate support beam.


Fig. 351: Cable replacement

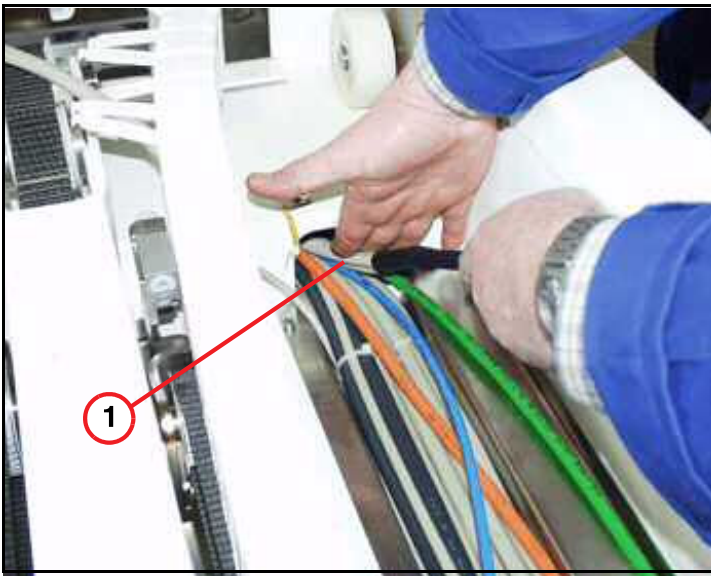

(1) Push the new cable

- At the same time as you pull the cable from the magnet side: Push the new cable from the table side through the tabletop plate support beam.

Fig. 352: Replace defective cable

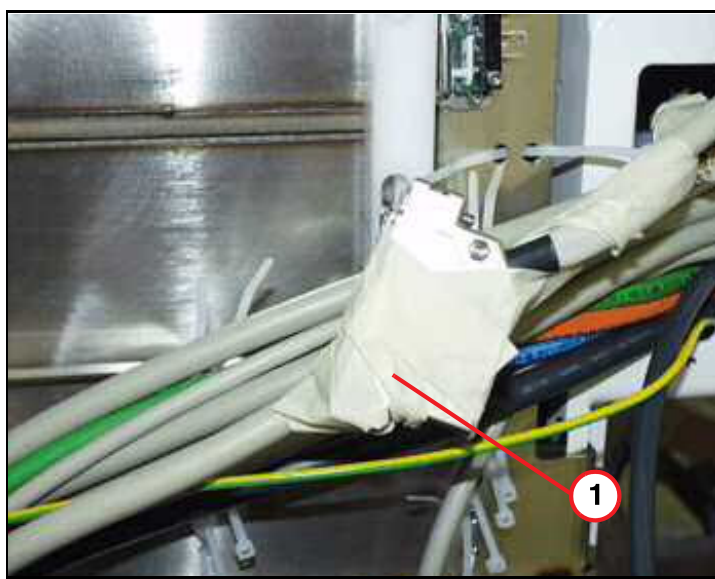

(1) Pull the new cable by using the defective one

- Feed the cable through the support beam. Pull the cable out until it is of the same length as the old one. Measured from end to end, that is, from the end of the energy chain to the end of the affected cable at table side.

Fig. 353: Energy chain for lifting column



- The replaced cable has to run parallel without crossing any other cable in the energy chain.
- Run the cable through the energy chain on the lifting column. Please ensure that the cables in the energy chain run in parallel.

- Align the new cable and install it with tie wraps.
- Connect the cable to the affected component.
- Close the chain links of the energy chains.

3.14.4.5 Adjustments

For adjustment of the toothed belt, please refer to (→ Adjustments / Page 170).

3.14.4.6 Concluding steps

- Reinstall mechanical components.
- Reinstall the covers of the patient table.

- Install the corresponding covers.

3.14.4.7 Final checks

- Switch ON F24 (RFIS).
- Switch ON the corresponding circuit breaker.
- Do a check of the proper table movement, without any grinding or noise.
- Do a check of the proper function of the replaced cable.

3.15 Control panel

3.15.1 General

This section describes the replacement of the control panels, the display, and the keypad foil.

3.15.2 Tools and time required

3.15.2.1 Time

Approx. 20 min./ 25 min.

3.15.2.2 Tools

Non-magnetic tool kit

3.15.3 Preliminary

- Switch OFF F22 (RF room).

3.15.4 Disassembly


3.15.4.1 Control panel left and rear side

There are two control panels on the left and right front side of the magnet. As an option, a third control panel can be placed at the back of the magnet.

This section describes how to disassemble the control panel on the left front side as well as the one in back.

Avanto


Fig. 354: Locking the 3. control unit into place

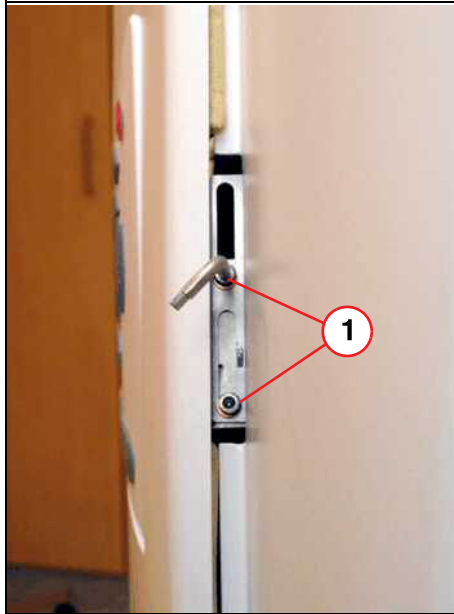

(1) Screws

- Remove ring segment of cover.
- Loosen both Allen screws, do not remove them.
- Shift the Allen screws up and carefully slide out the control panel.

Fig. 355: Installing the panel

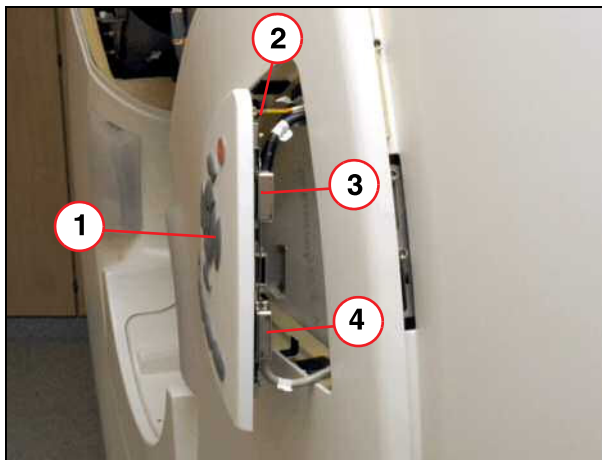

(1) 3. control unit
(2) PE (protective ground)
(3) X80
(4) X82

- Remove the control panel cables.

Espreo

Fig. 356: Control panel of Espreo



- To remove the left control panel, remove the screws on top of the service access cover.

- Remove the service access cover.

Fig. 357: Service access cover Espree



- Remove the screws of the defective control panel in back or the left front of the magnet.

Fig. 358: Mounting of left control panel



(1) Allen screws to be removed

- Carefully slide out the control panel.

Fig. 359: Disassemble control panel



- Disconnect the cables and remove the control panel.

Fig. 360: Disconnect cables of control panel




3.15.4.2 Control panel, right side

This section describes how to disassemble the control panel on the right side.

- Remove the lower right cover.

Fig. 361: Lower right cover

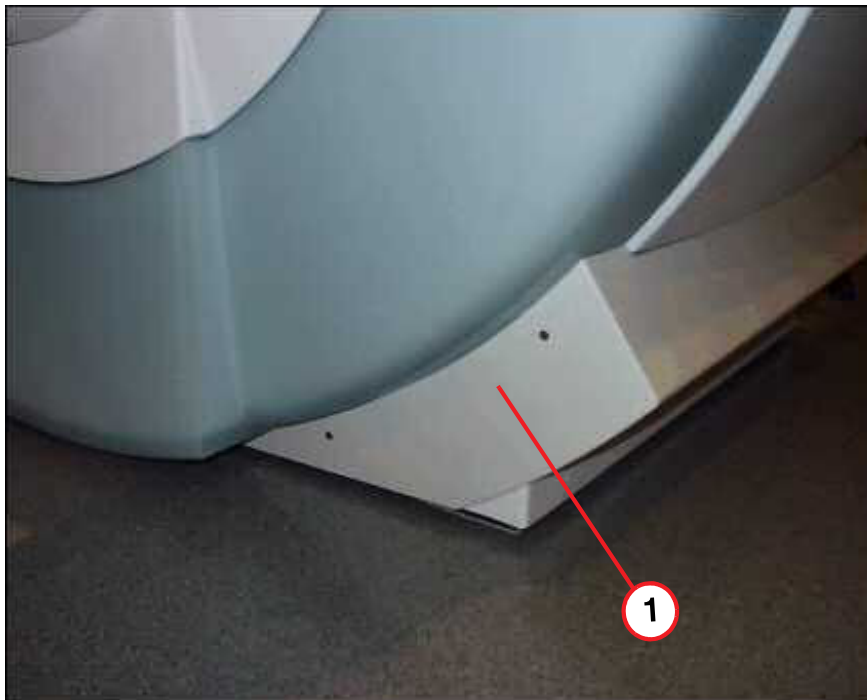

(1) Lower right cover to be removed

- Remove front ring segment to access the Allen screws at the right front display.
- Remove the Allen screws at the right front display.
- Carefully slide out the control panel.

Fig. 362: Remove Allen screws



- Disconnect the cables and remove the control panel.

Fig. 363: Disconnect cables of control panel



3.15.4.3 Keypad foil / Display electronics

This section describes how to disassemble the control panel. It is possible to replace the control panel electronics or the keypad foil.

- Remove the control panel and place it on a clean workbench.

- Remove countersunk screws and remove the display electronics.

Fig. 364: Display electronic removed



- If necessary, replace the display electronics or the keypad foil.

Fig. 365: Display electronic of control panel



Fig. 366: Keypad foil of control panel



3.15.5 Assembly

- Replace defective part.
- Reassemble the control panel.
- Reconnect the cables.

3.15.6 Concluding steps

- Remount the control panel into the cover
- Reassemble the covers and ring segment if necessary

3.15.7 Adjustments

- n.a.

3.15.8 Final checks

- Switch ON F22 (RF room)
- Test function of control panel

3.16 Patient light Espree

3.16.1 General

This section describes the replacement of the patient bore light located in the bore tube. There are two light rails in the funnel of the magnet (bore tube); they can be replaced separately.

3.16.2 Tools and time required

3.16.2.1 Time

Approx. 30 min.

3.16.2.2 Tools

Non-magnetic tool kit

3.16.3 Preliminary


- Switch OFF F22 (RF-room).

3.16.4 Disassembly

3.16.4.1 Funnel

Deinstalling the funnel


Fig. 367: Funnel

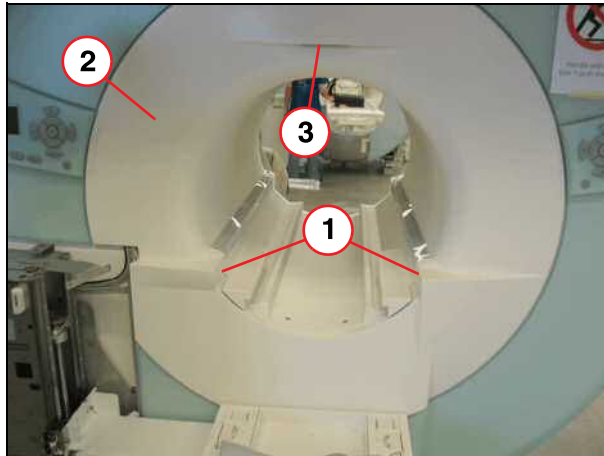

- (1) Screws
- (2) Funnel
- (3) Laser localizer faceplate

- Lower the patient table
- Remove the light marker faceplate (3).
- Remove the screws (1).

Fig. 368: Funnel, left side

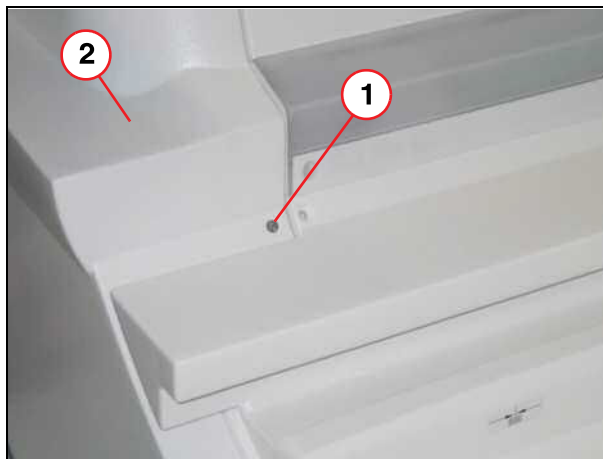

- (1) Screw
- (2) Funnel

Fig. 369: Remove funnel



(1) Funnel

- Remove the funnel (attached via Velcro).


Installing the funnel

- Attaching the covers, deinstallation in reverse order

Patient bore light

- Disconnect the cable from the defective patient bore light located at the front of the magnet.

Fig. 370: Cable for patient bore light

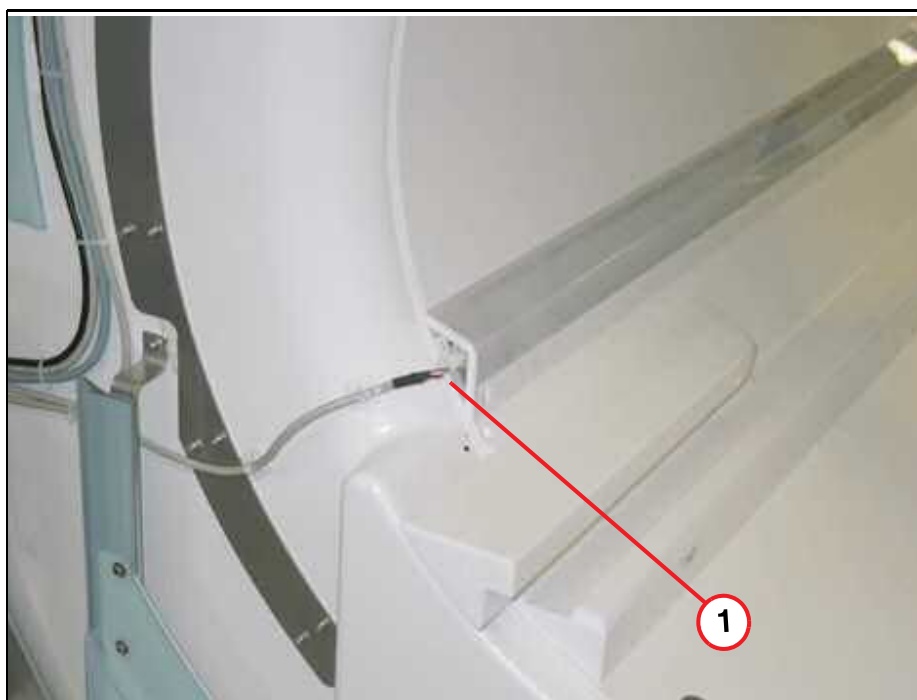

(1) Disconnect the cable

- Remove the screws of the defective light rail.

- Disconnect the cables in back of the light rail and remove them.

Fig. 371: Patient light in the funnel



(1) Screws to be removed

3.16.5 Assembly

- Replace defective part.
- Reconnect the cables.
- Reassemble patient bore light.

3.16.6 Concluding steps

- Reassemble the funnel in reverse order.
- Reassemble covers if necessary.

3.16.7 Adjustments

- n.a.

3.16.8 Final checks

- Switch ON F22 (RF-room).
- Test function of patient bore light.

3.17 Patient light Avanto

3.17.1 General

This section describes the replacement of the patient light located on the rear cover.

3.17.2 Tools and time required

3.17.2.1 Time

Approx. 10 min.

3.17.2.2 Tools


Non-magnetic tool kit

3.17.2.3 Preliminary

- Switch OFF F22 (RF-room).

3.17.3 Disassembly


Fig. 372: Rear cover

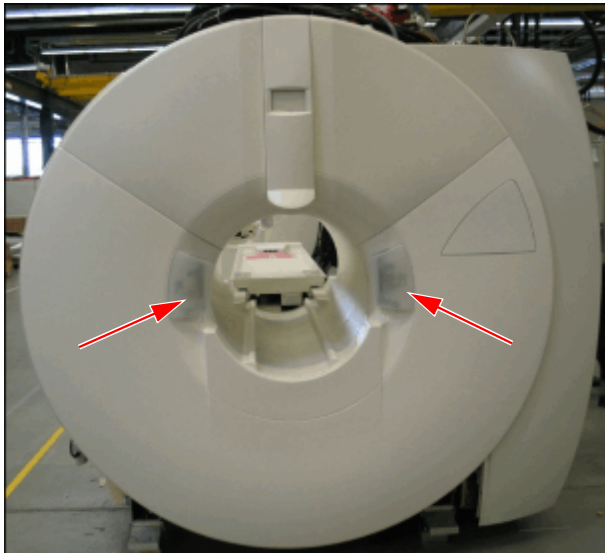

(1) LED cover

- Remove the screw on the upper side on the LED cover.

Fig. 373: Tunnel light

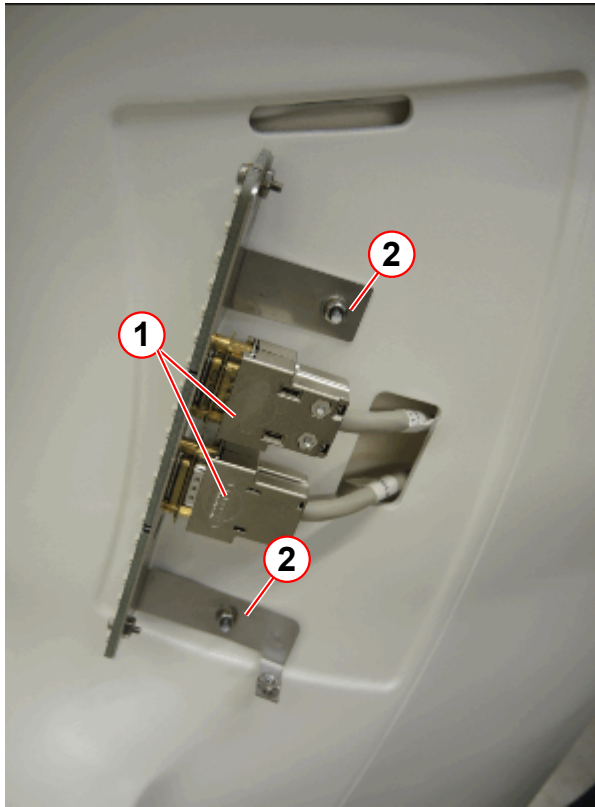

- (1) Plugs
(2) Mounting plate

- Unscrew and unplug the plugs (1).
- Loosen the nuts (2).

3.17.4 Assembly

- Replace defective part.
- Reconnect the cables.

3.17.5 Concluding steps

- Reassemble the LED cover.

3.17.6 Adjustments

- n.a.

3.17.7 Final checks

- Switch ON F22 (RF-room).
- Test function of patient bore light.

3.18 Light localizer / PMU antenna / microphone

3.18.1 General

This section describes the replacement of the components (light localizer, PMU antenna, microphone) located behind the front funnel.

There will be two microphones for future **Espreo** systems; they will be located at the front and the back of the magnet bore.

3.18.2 Tools and time required

3.18.2.1 Time

Approx. 30 min.

3.18.2.2 Tools

Non-magnetic tool kit

3.18.3 Preliminary


- Switch OFF F22 (RF-room).

3.18.4 Disassembly

3.18.4.1 Funnel

Deinstalling the funnel


Fig. 374: Funnel

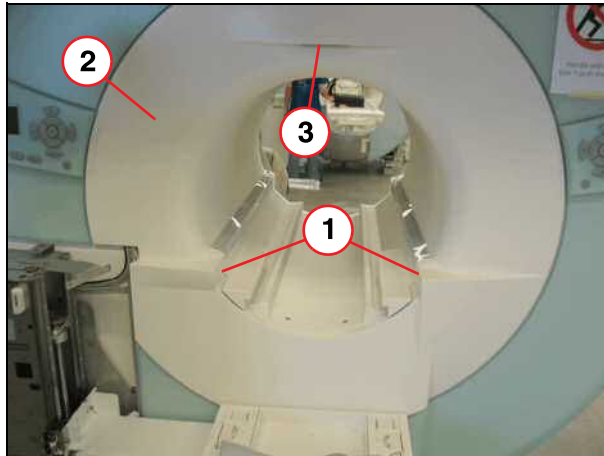

- (1) Screws
- (2) Funnel
- (3) Laser localizer faceplate

- Lower the patient table
- Remove the light marker faceplate (3).
- Remove the screws (1).

Fig. 375: Funnel, left side

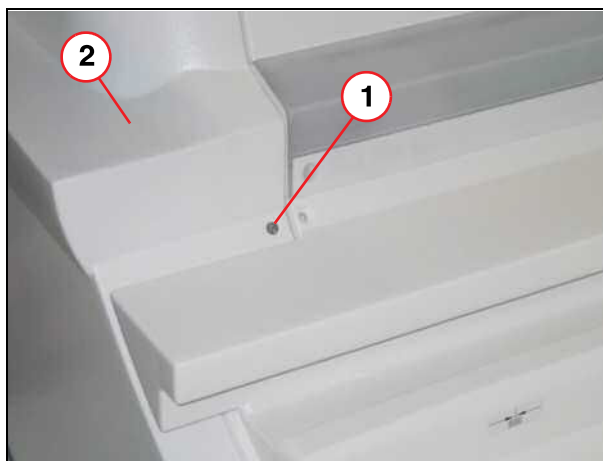

- (1) Screw
- (2) Funnel

Fig. 376: Remove funnel



(1) Funnel

- Remove the funnel (attached via Velcro).

Installing the funnel


- Attaching the covers, deinstallation in reverse order

3.18.4.2 Light localizer, microphone, PMU antenna

This section describes how to disassembly the components located on the assembly plate behind the front funnel. It may also be necessary to have access to these components for troubleshooting.

- If you have to remove the assembly plate, i.e. for troubleshooting, mark the position beforehand.

Fig. 377: Assembly plate

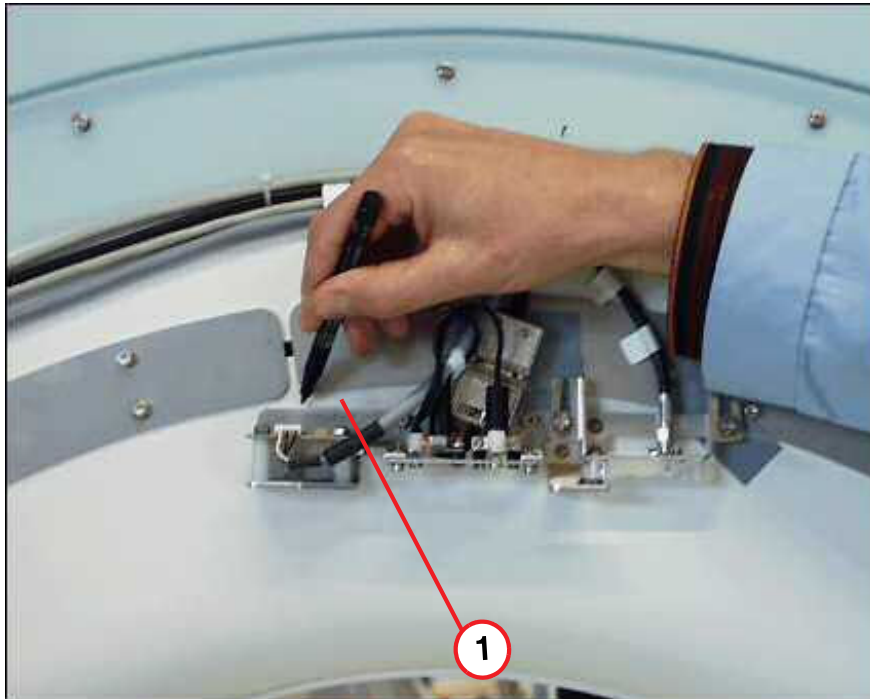

(1) Mark the position of the assembly plate

- If necessary cut the tie wraps, disconnect, and remove cables.

Fig. 378: Assembly plate, light localizer



- Disassembly the defective part.

3.18.5 Assembly


- Replace defective part.
- If necessary, reroute cables along their original routing.
- Reassemble the assembly plate.
- Reconnect the cables.

3.18.6 Concluding steps

- Reassemble the funnel in reverse order.
- Reassemble the covers.

3.18.7 Adjustments

Fig. 379: Light localizer

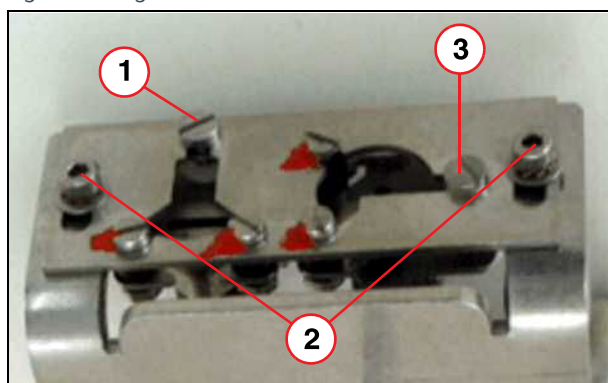

- (1) Setting screw A
- (2) Attachment B
- (3) Setting screw C


The principle for adjusting the light localizer is the same for both the Avanto and the Es-
pree.

Adjustment screw A (1), for laser adjustment A horizontal .

Mounting B (2), for laser adjustment B horizontal.

Adjustment screw C (3), for laser adjustment B vertical.


Fig. 380: Laser 1

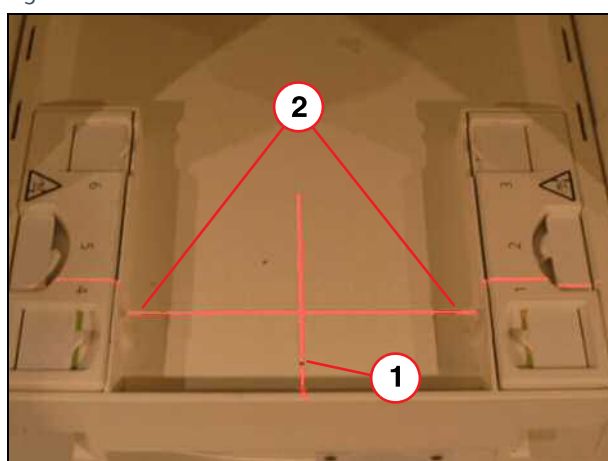

- (1) Laser A horizontal
- (2) Laser B horizontal

- Release the two screws and adjust laser B (2) horizontally to the tabletop plate, the left and right distance are equal.
- Use screw A and adjust laser A (1) horizontally to the reference hole in the tabletop plate.

Fig. 381: Laser and auxiliary tools

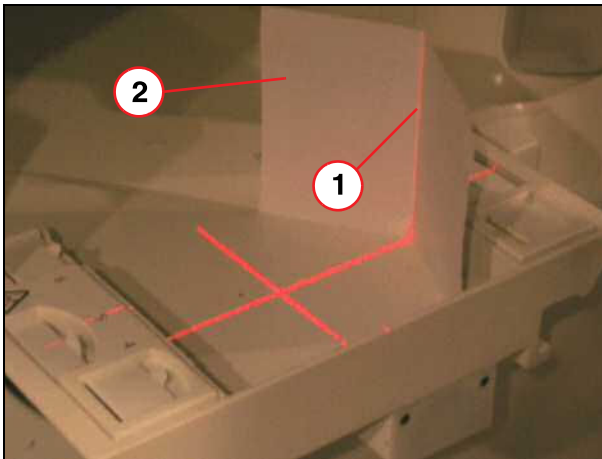

- (1) Laser B vertical
- (2) Auxiliary tool, folded paper

- Crumple up a white piece of paper and place it on the tabletop plate.
- Use screw C and adjust laser B (2) vertically along the white paper.

3.18.7.1 Concluding steps

- After finishing the adjustment, remount the laser cover.

3.18.8 Final checks

- Switch ON F22 (RF-room).
- Do a test of the function of the replaced part.

4.1 Locking Hooks (Mechanical)

4.1.1 General

The Trolley consists of several components. Due to the mechanical complexity of the trolley, the replacement of FRUs is limited to locking mechanics, PCBs and keypads described in this chapter. In case the hydraulic components need to be replaced, the trolley has to be returned to the factory.

This chapter describes how to replace locking hooks. In addition, this procedure replaces several other mechanical procedures, e.g. raster element, spring element, Bowden wire....

4.1.2 Tools and Time required

4.1.2.1 Time

Approx. 120 min

4.1.2.2 Tools

Non-magnetic tool kit


4.1.3 Preliminary

- Move the trolley out of the examination room.

4.1.4 Disassembly

4.1.4.1 Covers


Fig. 382: Trolley, front cover

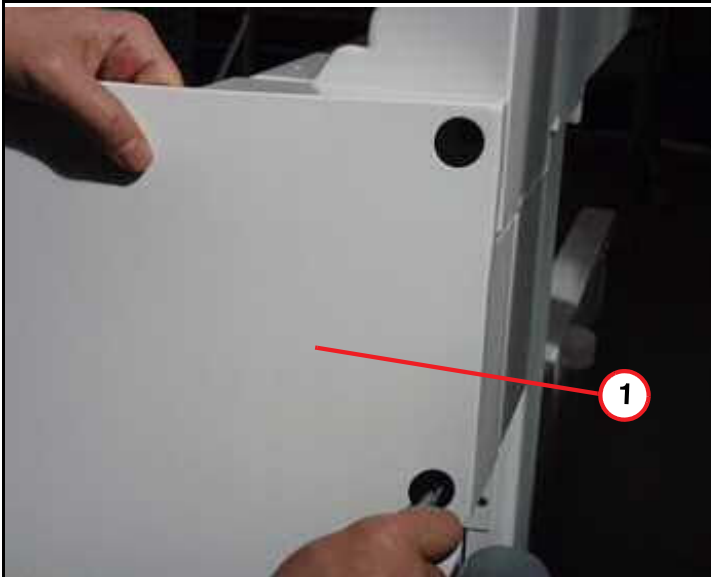

(1) Front cover to be removed

- Remove the 4 plastic caps and remove the front trolley front cover by taking off the four Allen screws.


Fig. 383: Lower cover

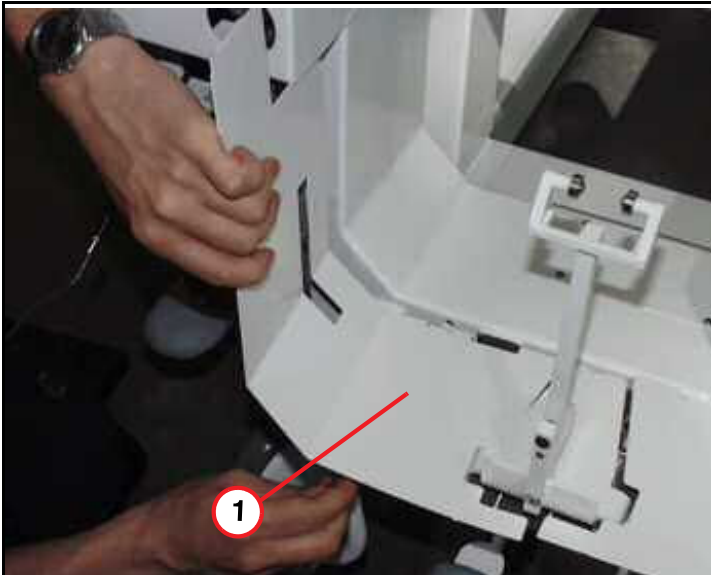

(1) Lower cover to be removed

- Remove lower cover by taking off the Allen screws.


Fig. 384: Allen screw of front cover

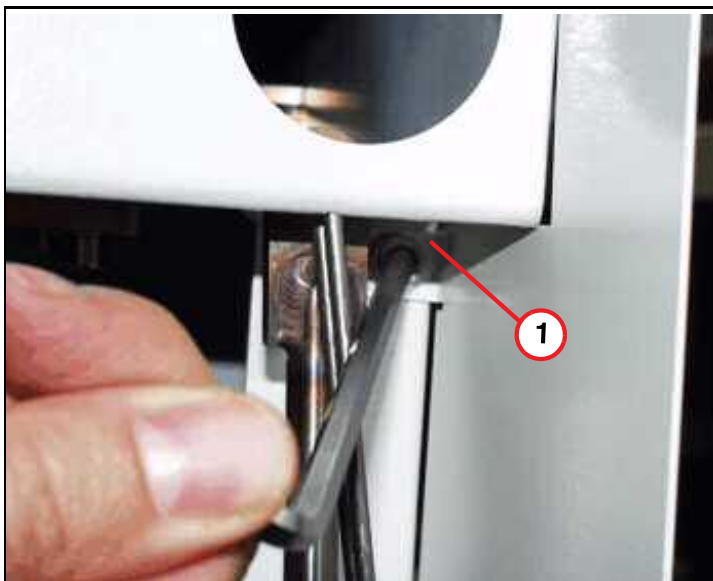

(1) Allen screw to be removed

4.1.4.2 Locking mechanics


Fig. 385: Clamping plate

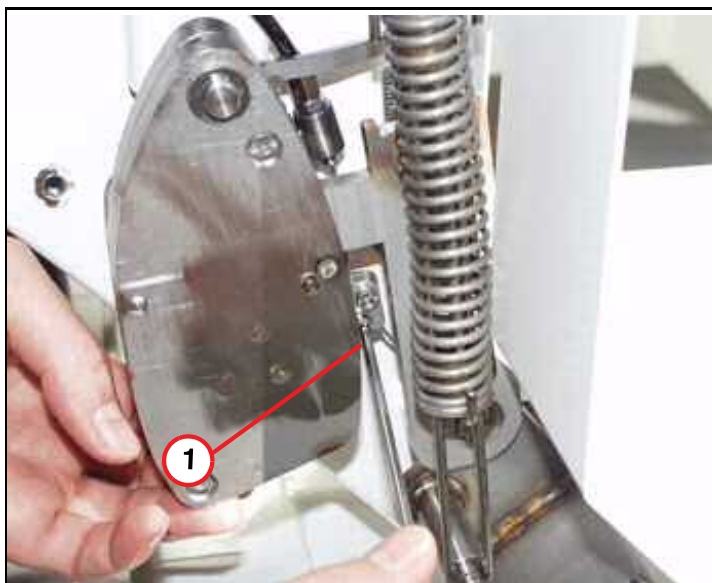

(1) Clamping plate of Bowden wire

- Release clamping plate from the Bowden wire and remove the wire.


Fig. 386: Locking element

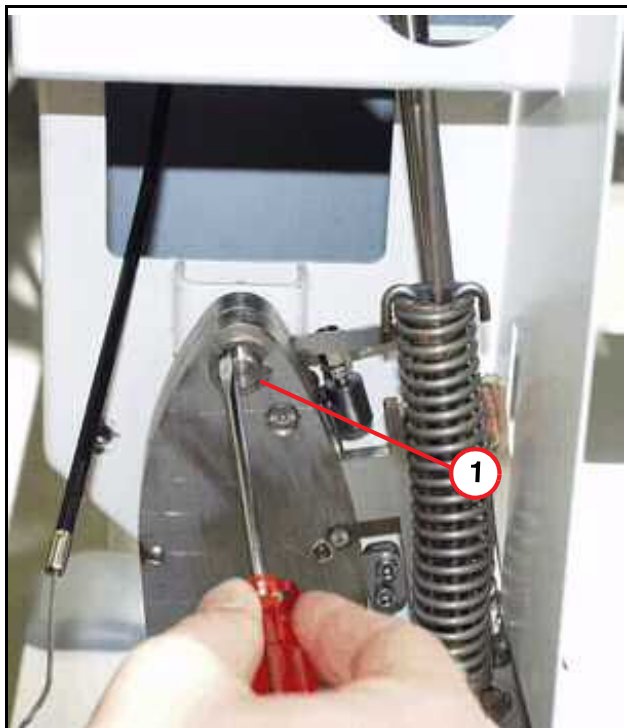

(1) Split ring of locking element


- Remove split ring of rast element.

Fig. 387: Rast element

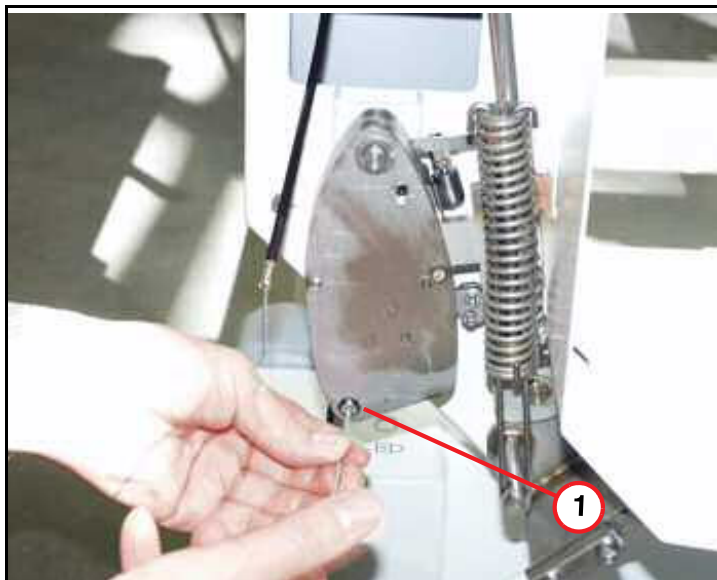


(1) Allen screws of rast element

- Remove the two Allen screws.

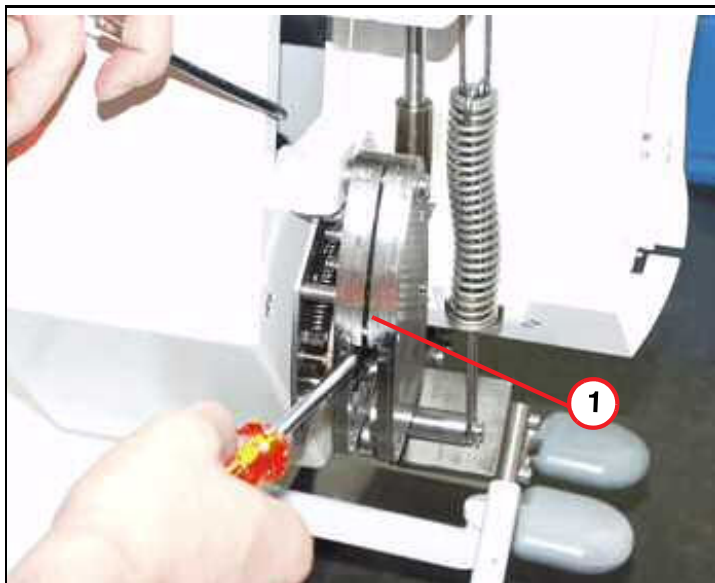Fig. 388: Rast element

(1) Rast element

- Remove the front plate of the lock-in element by using a flat head screwdriver.

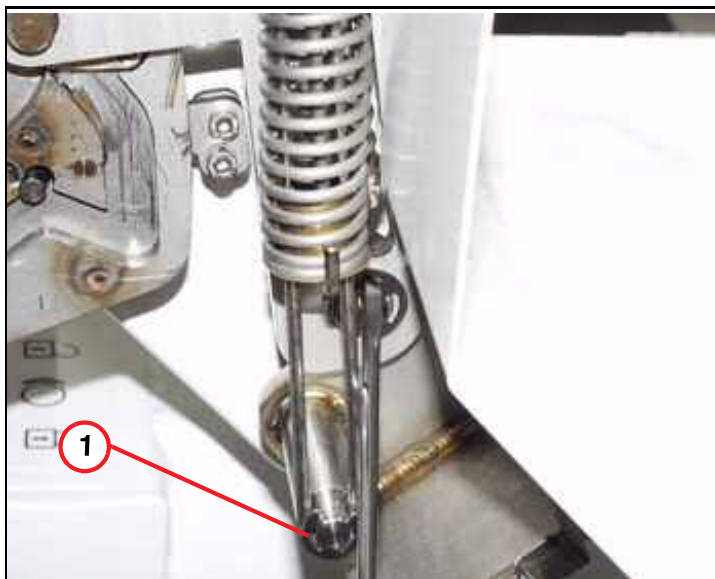
Fig. 389: Spring element

(1) Lower split ring of spring element

- Remove the split ring.

Fig. 390: Actuating spring


- Remove bolt and spring element.

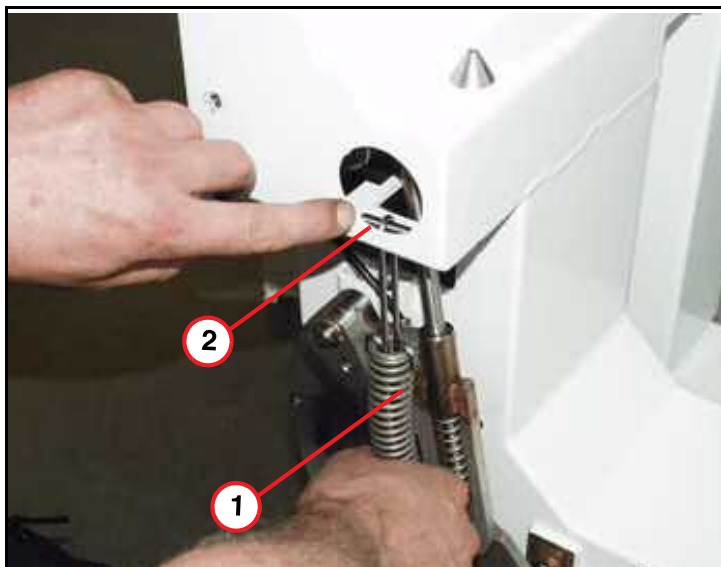

- (1) Spring element
- (2) Bolt to be removed


Fig. 391: Split ring

- Remove upper split ring.

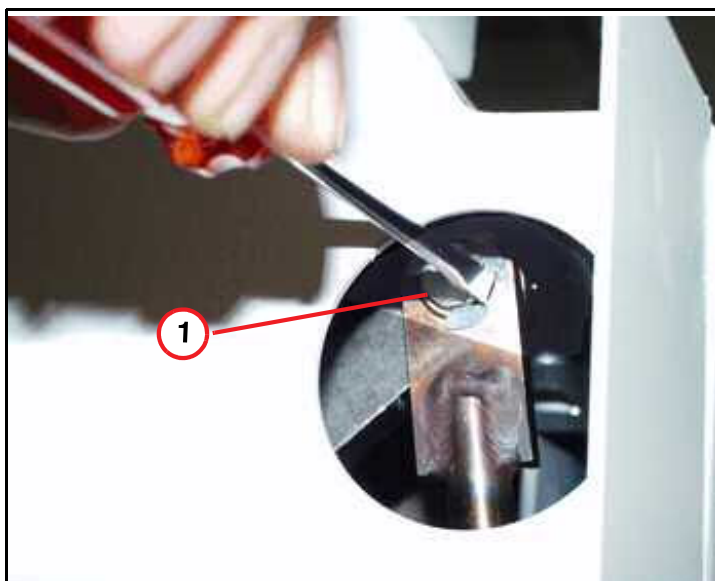

- (1) Upper split ring

Fig. 392: Spring rod

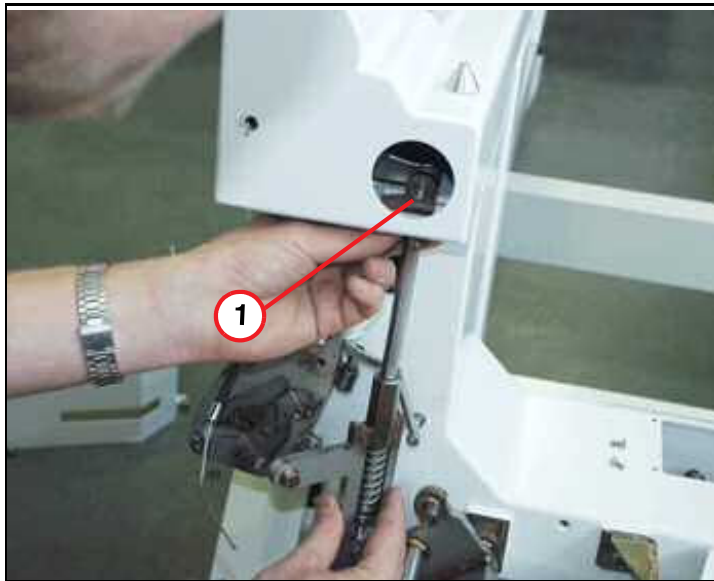

(1) Spring rod

- Remove lower split ring and the spring rod.


Fig. 393: Grub screw



(1) Grub screw of actuating lever

- Remove the grub screw of the actuating lever.


Fig. 394: Actuating lever

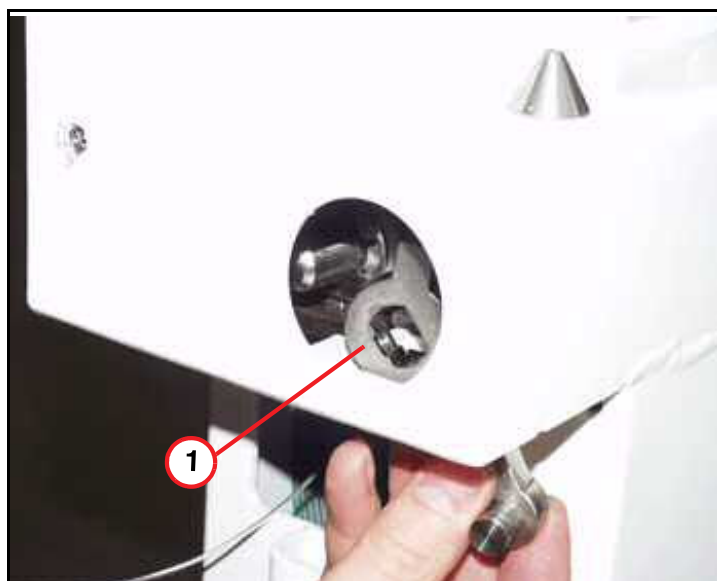

(1) Actuating lever removed

- Remove actuating lever.

Fig. 395: Trolley locking mechanics

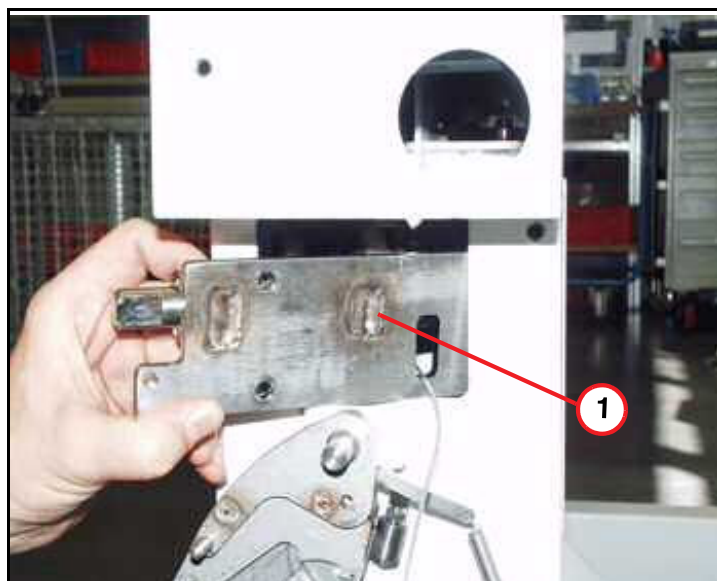

(1) Trolley locking mechanics removed

- Remove the four Allen screws.
- Remove locking mechanics and replace them.

4.1.5 Assembly

4.1.5.1 Locking mechanics

Fig. 396: Spring rod



(1) Spring rod

- Remount the actuating lever and tighten the grub screw.
- Remount spring rod.


Fig. 397: Bowden wire for tabletop locking



(1) Insert Bowden wire bush

- Reassemble the bushing of the Bowden wire and do a check of the function of the locking mechanism.

Fig. 398: Clamping plate

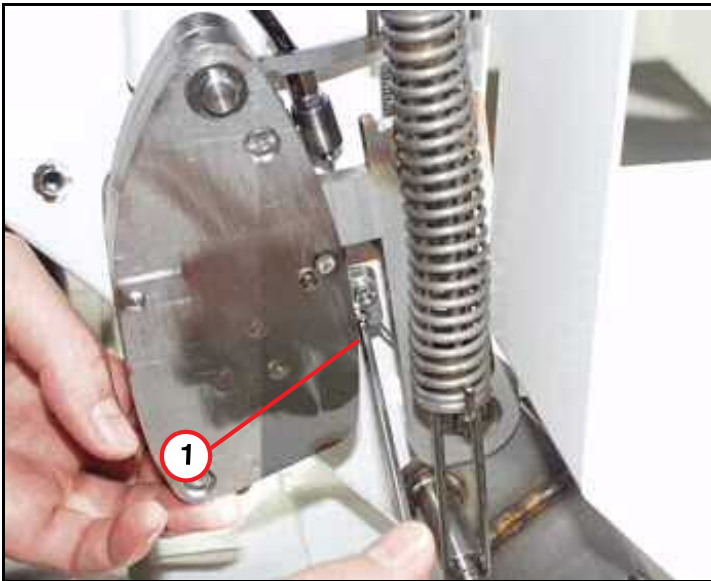

(1) Clamping plate of Bowden wire

- Insert the Bowden wire into the clamping plate, slowly pull Bowden wire until the plastic hook reacts. Release it a little bit and remain in this position, tighten the clamping plate.


Fig. 399: Trolley hook



(1) Trolley hook assembly

- Insert the split rings from the spring rod
- A properly adjusted Bowden wire leads to the planned hook assembly (plastic and metal hook are level). Tighten the counter nut after the adjustment.

Fig. 400: Rast element

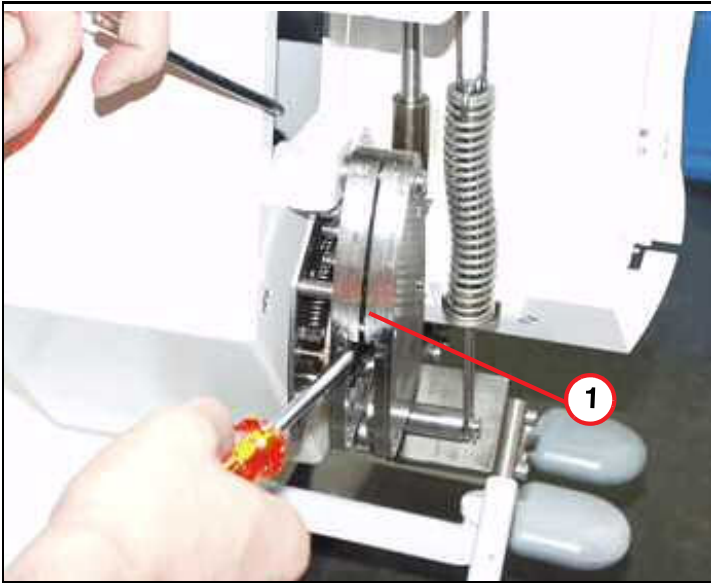

(1) Rast element

- Reassemble the front plate of the lock-in element and insert the Allen screws, the standoff and the split ring.
- Install the spring element, with bolt and split ring.

4.1.6 Final Steps

- Reinstall the covers.

4.1.7 Final Checks

- Position the trolley at the patient table.
- Start trolley function and check for smooth and safe operation.

4.2 Keypads / PCBs

4.2.1 General

This chapter describes how to replace locking hooks. In addition, this procedure replaces several other mechanical procedures, e.g. lock-in element, spring element, Bowden wire....

4.2.2 Tools and time required

4.2.2.1 Time

Approx. 20 min.

4.2.2.2 Tools

Non-magnetic tool kit

4.2.3 Preliminary

- Move the trolley out of the examination room.

4.2.4 Disassembly

4.2.4.1 Keypads / PCBs

Fig. 401: Cover of keypad



- Remove the small cover of the left or right side.


(1) Cover to be removed

Fig. 402: Battery cover, trolley remote control



- Remove the cover of the PCB.


Fig. 403: Circuit board

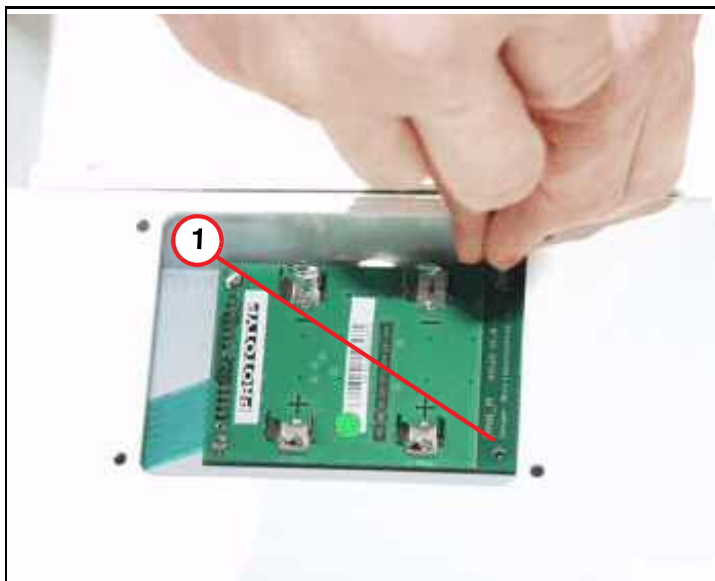

- Remove batteries and PCB.

(1) Screws of PCB removed


Fig. 404: Cable connections

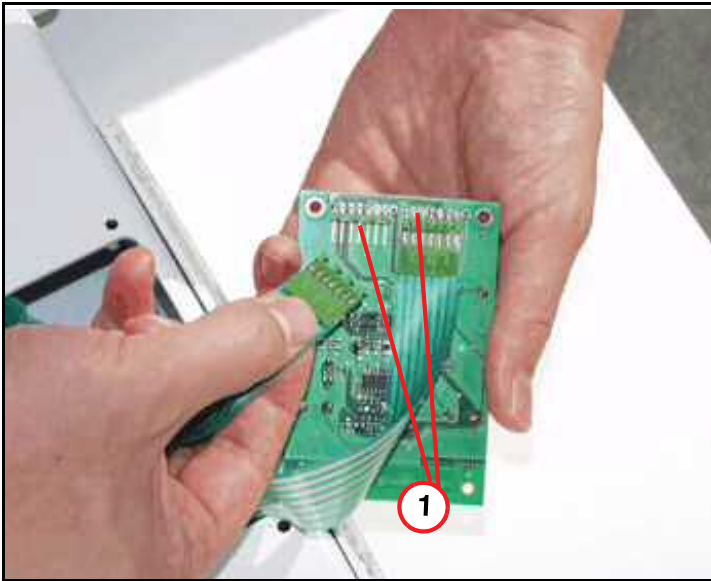

(1) Cable connections on PCB

- Disconnect cables and replace PCB, if necessary.

Fig. 405: Keypad

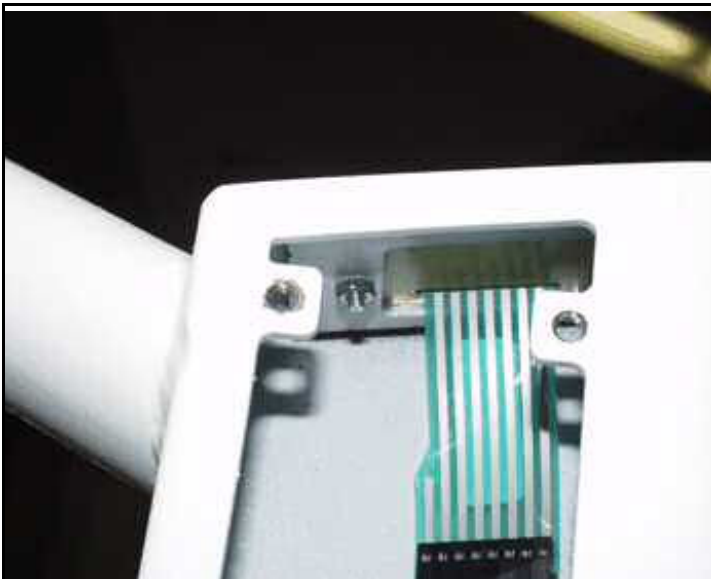

- Replace the defective keypad.

4.2.5 Assembly

- Reassemble in reverse order.

4.2.6 Concluding steps

- Reinstall the covers.

4.2.7 Final checks

- Position the trolley at the patient table.
- Start trolley function and check remote control buttons

4.3 Wheel Reinforcement Kit

4.3.1 General

This chapter describes how to install the reinforcement kit (material number 10092160). By installing a wheel reinforcement set in the trolley, the mechanical stress resistance of the trolley wheels improves.

4.3.2 Tools and Time required

4.3.2.1 Time

60 min

2 Persons required

4.3.2.2 Tools


- Heat gun
- Set of Allen keys
- Flat wrench 13 mm
- Loctite 577 (material number 03443348)
- Safety gloves
- Permanent waterproof pen
- Isolation tape
- Protect the trolley surface when putting it on its side (e.g. with a blanket, weight of the trolley is approx. 150 kg).

4.3.3 Preliminary

- Move the trolley out of the examination room.

4.3.4 Installing the Wheel Reinforcement Kit


Fig. 406: Wheel Reinforcement Kit

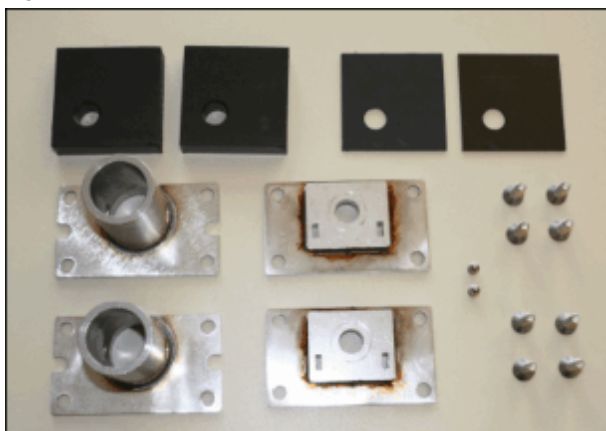

- Check the scope of delivery.

4.3.4.1 Rear Wheels.


Fig. 407: Removing end stop

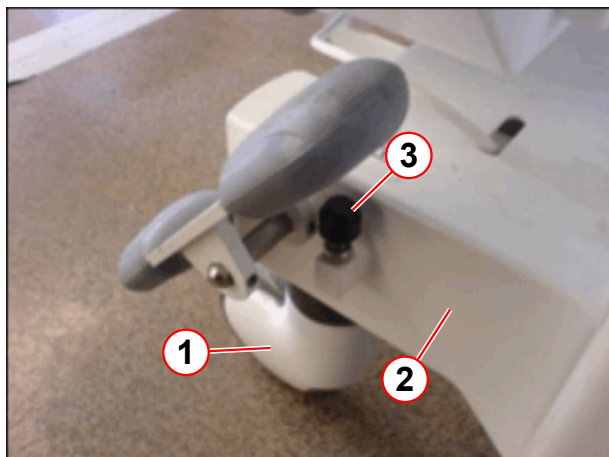

- (1) Rear wheel
- (2) Base frame
- (3) End stop

- Move the trolley forward (approx. 0.5 m) so that the rear wheels (1) are parallel to the chassis (2).
- Activate the brakes.
- Loosen the counter nut at the end stop (3)
- Remove the end stop (3).


Fig. 408: Brake pedal

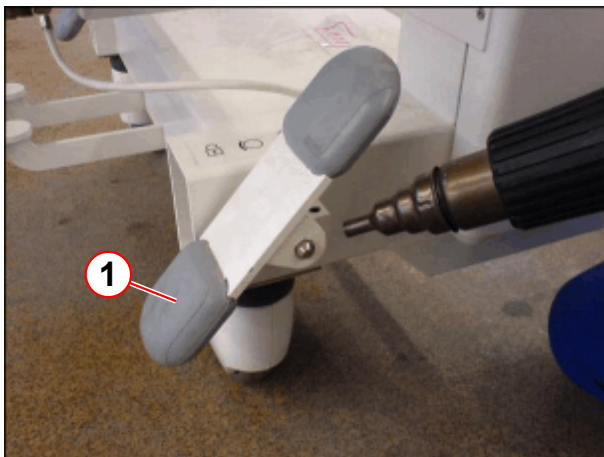

(1) Brake pedal

- Heat up the threaded bolts of the brake pedal (1) with a heat gun.
- Be careful not to damage the heat-sensitive surface.

Fig. 409: Brake

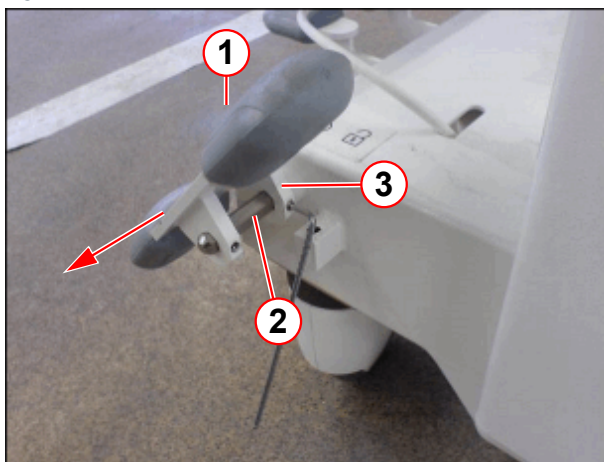

(1) Pedal
(2) Hexagonal rod
(3) Spacer

- Take an Allen key and remove the threaded bolts of the pedal (1)
- Pull out pedal (1) and spacer (3) (white plastic)




The rear wheels must be parallel with the base frame and the brake must still be engaged.

Fig. 410: Trolley on the side



- Put the trolley on the side, where the pedal has been removed (cover floor to protect the surface of the trolley, trolley weight approx. 150 kg).
- The following describes the replacement of the wheel take-ups for the front and rear wheels.


Fig. 411: Pulling out the hexagonal rod

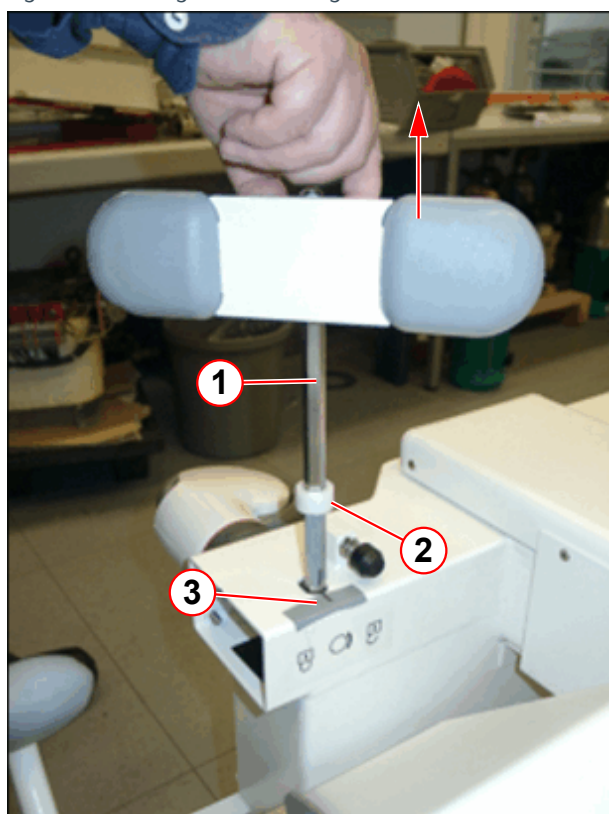

- Mark the position of the pedal angle at the rod (3) (e.g. on a strip of isolation tape) – for correct reinstallation.
- Pull the pedal out completely together with the hexagonal rod (1) and spacer (2)

- (1) Pedal with hexagonal rod
 (2) Spacer
 (3) Marking the pedal angle position


Fig. 412: Screws on the rear wheel

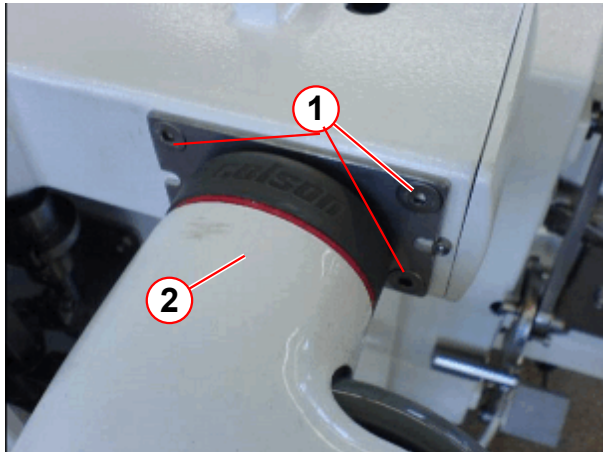

- (1) Screws
(2) Wheel with wheel take-up

- If necessary, carefully heat up the screws (1) of a rear wheel
- Take off the rear wheel (2) with wheel take-up


Fig. 413: Wheel take-up

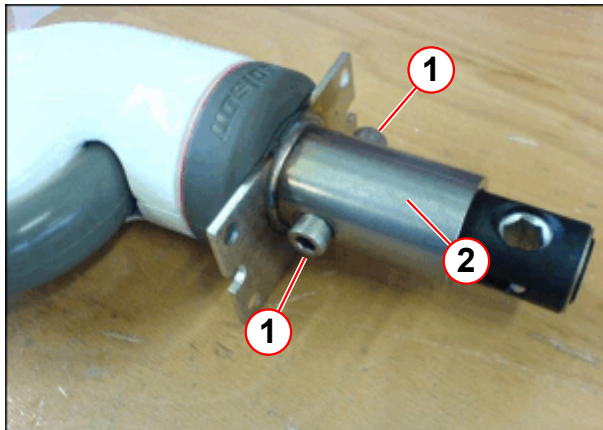

- (1) Screws
(2) Wheel take-up

- The wheel take-up is not symmetrical.
- Fix the wheel take-up (2) with the 2 screws. (1)


Fig. 414: Reinstall the wheel

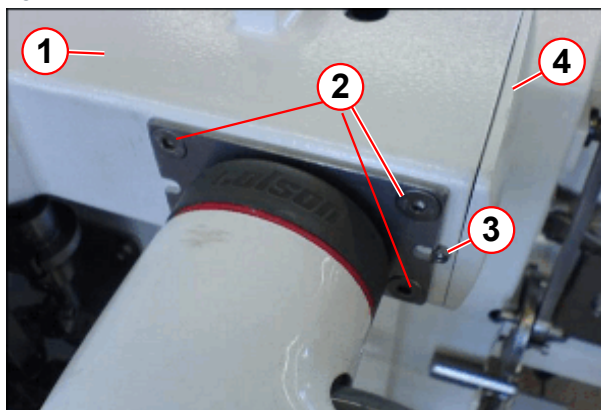

- (1) Base frame
- (2) Screws (4)
- (3) Screw
- (4) Cover

- Install the modified rear wheel on the base frame (1) with the leveling plates.
- Apply Loctite 577, P.N. 03443348, to screws (2) and tighten them
- Loosen the screw (3) for the cover (4) approx. 8 mm
- Repeat the procedure for the other rear wheel

Fig. 415: Taking off the cover

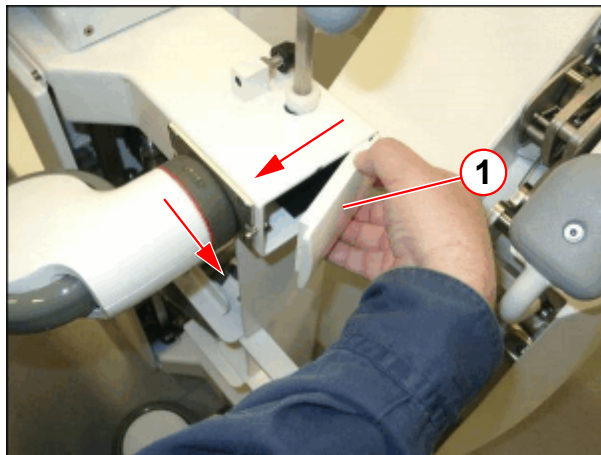


- (1) Cover

Remove the cover (1) as follows:

- Shift the cover toward the base frame
- Pull out cover


Fig. 416: Installing distance plate

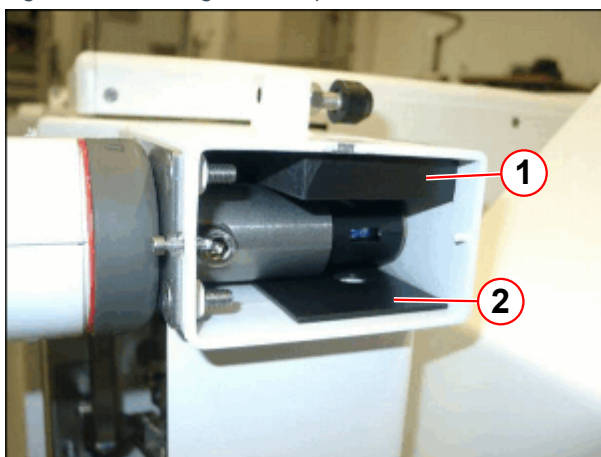

- (1) Distance plate 15 mm
- (2) Distance plate 2 mm

- Install distance plates, 15 mm (1) and 2 mm (2).
- Repeat the procedure for the other rear wheel.

Fig. 417: Pulling in the hexagonal rod

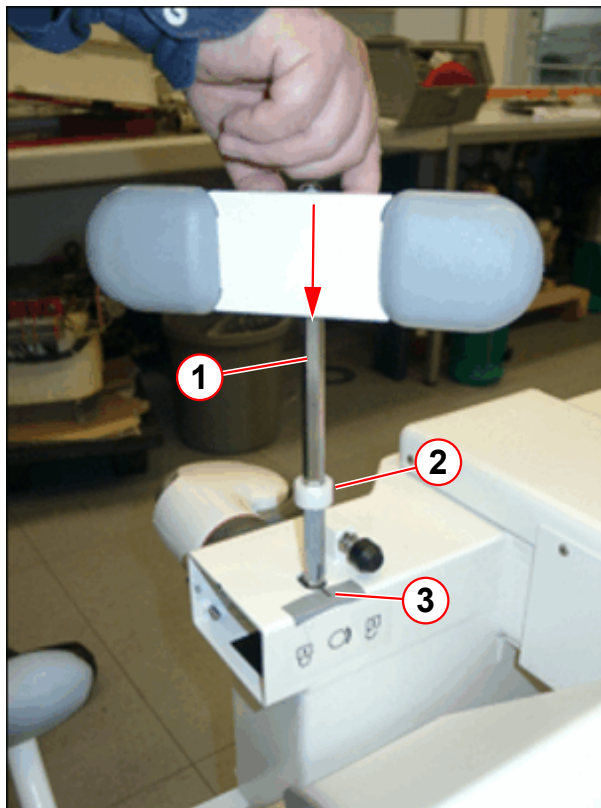


- (1) Pedal with hexagonal rod
- (2) Spacer
- (3) Marking the pedal angle position

- Slide the hexagonal rod with pedal (1) and spacer (2) back in. Make sure that the rod is sliding through the holes of the distance plates.
- After correctly inserting the rod, remove marks (3)
- Install the cover and secure it with the screw

4.3.4.2 Front Wheels


Fig. 418: Remove the front wheels

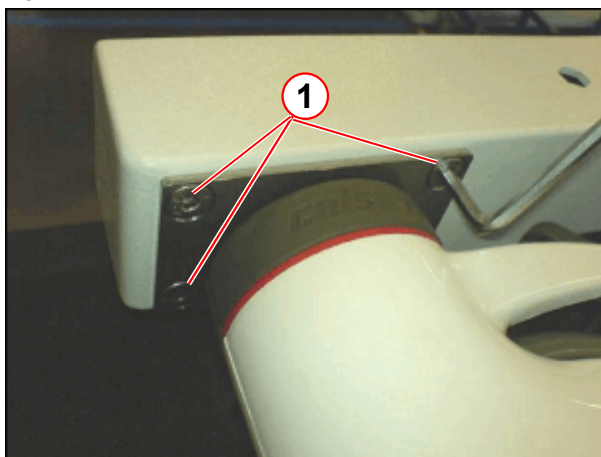

- (1) Screws

- If necessary, carefully heat up the screws (1) of both front wheels
- Take off the front wheels

Fig. 419: Replace wheel mount

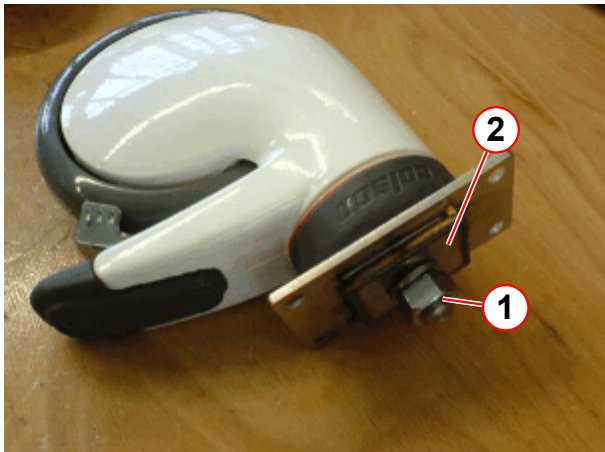

- (1) Nut
- (2) Mount


- Open the nut (1) and remove the wheel mount.
- Install the new wheel mount in the same way as the old one (material is thicker); use Loctite 577.
- Tighten the nut (1) securely.
- Repeat the procedure for the other front wheel.

Fig. 420: Front wheel



- Tighten/loosen the screw on the front wheel plates shown in the picture.


Fig. 421: Install the front wheel

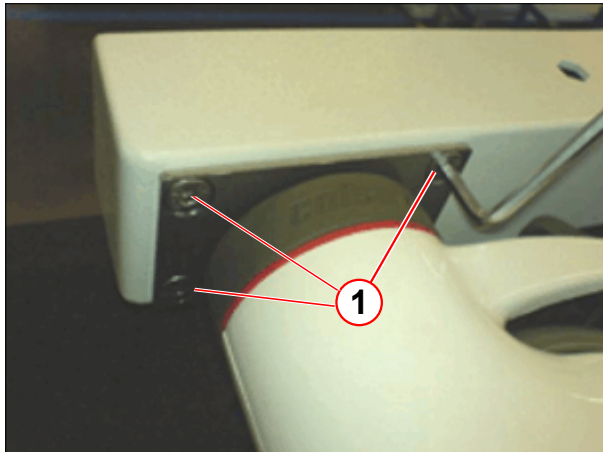

(1) Screws


- Install the modified front wheel with the leveling plates.
- Tighten screws (1) securely, use Loctite 577.

Fig. 422: Lift up trolley

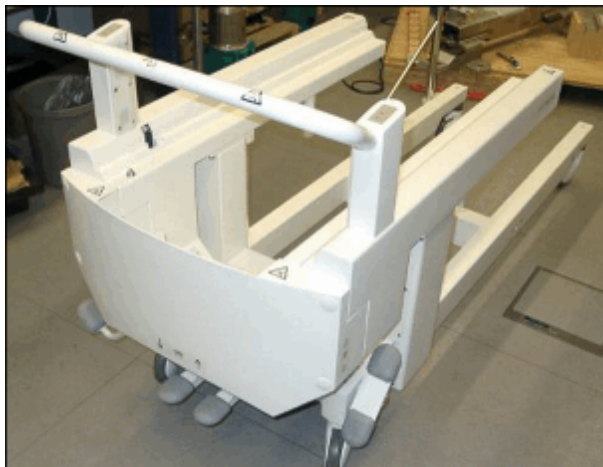

- Lift the trolley carefully. Be sure that the brake is still active.


Fig. 423: Installing pedal and end stop

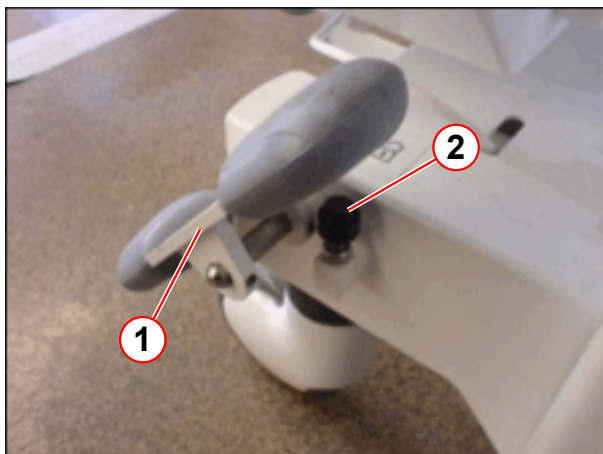

(1) Pedal
(2) End stop of brake pedal

- Reinstall the pedal (1) using Loctite 577 for the threaded bolts.
- The 2 pedals must be parallel
- Insert and align the end stop (2) of the brake pedal. Both pedals must be horizontal when not braked.

Fig. 424: Revision level

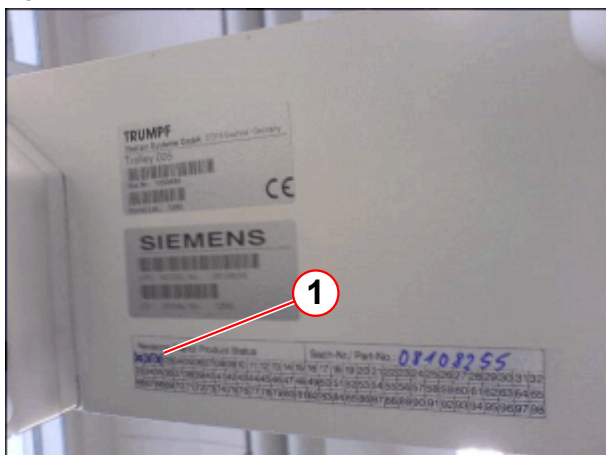

(1) Revision level label

- Mark the new revision level (refer to table below) on label (1) using a permanent, waterproof pen

Trolley System designation	Material number	Revision level
Trolley 005	08108255	03
Trolley 094	10018482	03
Trolley 092	10018397	02
Trolley 005/ 022 - 100 mm increased	10018902	01

4.3.5 Final Checks

- Position the trolley near the patient table.
- Start trolley function and check remote control buttons.

5.1 PDAU

5.1.1 General

The PMU consists of several components. However, only the exchange of PDAU needs to be described.



With software version VB13, a new PDAU has been introduced.

Older software versions will not support the new PDAU, material number 10096522, refer to spare part list.

5.1.2 Tools and time required

5.1.2.1 Time

Approx. 20 min.

5.1.2.2 Tools

Non-magnetic tool kit

Non-magnetic 2.5 mm Allen key

5.1.3 Preliminary steps

- Switch OFF F22 (RF-room).

5.1.4 Disassembly

5.1.4.1 Covers

- Remove the magnet side cover and all covers from the lifting column.
- Remove the cover for the patient table electronics.

5.1.4.2 PDAU

Fig. 425: PDAU, top of ptab-electronics



- Disconnect all cables from the PDAU.

- Remove the Allen screws and replace the PDAU.



The PDAU consists of small magnetic objects. When replacing the PDAU, observe caution. Maintain a safe distance to the magnetic field during this procedure.

Fig. 426: Remove Allen screws of PDAU



5.1.5 Assembly

- Insert the new PDAU and tighten the screws.

5.1.6 Concluding steps

- Reconnect the cables.
- Reassemble the covers.

5.1.7 Adjustments

- n.a.

5.1.8 Final checks

- Switch ON F22 (RF-room).
- Perform, reboot scanner.
- Check PDAU for proper function.

6.1 Loudspeaker/Microphone

6.1.1 General

This section describes the replacement of the loudspeaker and microphone on the Avanto.

6.1.2 Tools and time required

6.1.2.1 Time

Approx. 20 min. / 30 min.

6.1.2.2 Tools

Non-magnetic tool kit

6.1.3 Preliminary

- Switch OFF F22 (RF-room).

6.1.4 Disassembly

6.1.4.1 Covers/Loudspeaker

Fig. 427: Magnet side cover



- Remove the magnet side cover.

Fig. 428: Cover plate



- Loosen the screws and remove the cover plate.

Fig. 429: Foam rubber of loudspeaker



- Disconnect the cable and remove the foam rubber of the loudspeaker.

6.1.5 Assembly

- Replace loudspeaker.

6.1.6 Concluding steps

- Reconnect the cables.
- Reassemble the covers.

6.1.7 Adjustments

- n.a.

6.1.8 Final checks

- Switch ON F22 (RF-room).
- Test intercom unit.

6.1.9 Disassembly

6.1.9.1 Covers/Microphone

Fig. 430: Screws of small front funnel



- Move the patient table out of the magnet.
- Loosen (but do not remove) the 3 screws using a flat blade screwdriver.

Fig. 431: Packing ring of front cover



- Remove the packing ring of the front cover.

Fig. 432: Locking pin of small front funnel



- Pull the locking pin on the left and right front cover.

Fig. 433: Small front funnel



- Remove the small front funnel of the magnet front; pay attention to the microphone cable.
- Disconnect microphone cable.

6.1.10 Assembly

- Remove the microphone and replace it.

6.1.11 Concluding steps

- Reassemble the covers.

6.1.12 Adjustments

- n.a.

6.1.13 Final checks

- Switch ON F22 (RF-room).
- Test intercom unit.

6.2 CCD camera

6.2.1 General

This section describes the replacement of the CCD patient monitoring camera.

6.2.2 Tools and time required

6.2.2.1 Time

Approx. 15 min.

6.2.2.2 Tools

Non-magnetic tool kit

6.2.3 Preliminary

- Switch OFF F22 in the ACC.

6.2.4 Disassembly

6.2.4.1 Covers

Fig. 434: Screws for cover plate of camera



- Remove the two flat blade screws from the cover plate of the CCD camera.

Fig. 435: Screws for outer camera cover



- Remove the two Allen screws inside the outer camera cover.

Fig. 436: Back cover



- Carefully slide the camera cover up and remove them

6.2.5 Concluding steps

- Reconnect the cables.
- Reassemble the covers.

6.2.6 Adjustments

- n.a.

6.2.7 Final checks

- Switch ON F22 (RF-room).
- Visual test of patient monitoring.

- » The print number (M6-020.841.15) from the Replacement of Parts is revised with controlled language.
- » Editorial changes.

Version 07

- Hint: "Protective earth connection test" was added.
- Caution: "Hold the ferromagnetic part with both hands" was added.
- The assembly procedures for the partial drive MDSD and the Drive unit TK have been updated with a new mechanical key installation.
- Caution: "Possible damage of the spindle nut" was added.

Version 08

- Arresting cable added for exchanging and mounting the horizontal drive unit (→ Removal of horizontal drive unit / telescopic drive / Page 139)

Version 09

- Insert the wheel reinforcement kit procedure (→ Wheel Reinforcement Kit / Page 253)

Version 10

- Tools for guide rails: Allen key added (→ Tools / Page 16)
- Final test step, test spindle nut safety switch S30 added (→ Final checks / Page 33)
- Replacement tunnel light Avanto added (→ Patient light Avanto / Page 230)

Version 11

- Add new graphics (→ Spare part items of spindle / Page 39)
- Minor editorial changes (Charm MR 00447613)
- Add in the FRU: Transformer A4114 (→ FRU Overview / Page 10)

Version 12

- RF Coil Plugs (head end): (→ Assembly / Page 191)
Caution: "Any contact between the shield of the cable and the housing of the trap must be avoided."
and graphics (→ Fig. 317 Page 191), (→ Fig. 318 Page 192) and (→ Fig. 319 Page 193) added. (MR Charm 441285)

Version 13

- Minor editorial changes.

Version 14

- RF coil plugs: Required tool and service plug added.
- Brake disc: New photo added. Note added that the grooved side must face up. (→ Brake disk / Page 67)

There are no Hazard IDs in this document.

All documents may only be used for rendering services on Siemens Healthcare Products. Any document in electronic form may be printed once. Copy and distribution of electronic documents and hardcopies is prohibited. Offenders will be liable for damages. All other rights are reserved.

healthcare.siemens.com/services

Siemens Healthineers Headquarters
Siemens Healthcare GmbH
Henkestr. 127
91052 Erlangen
Germany
Telephone: +49 9131 84-0
siemens.com/healthineers

Print No.: M6-020.841.15.14.02 | Replaces: M6-020.841.15.13.02
Doc. Gen. Date: 05.17 | Language: English
© Siemens Healthcare GmbH, 2006

siemens.com/healthineers