

SIEMENS

MAGNETOM

MR

Software

System

Software Installation (VA25A)

© Siemens, 2004

All documents may only be used for rendering services on Siemens Healthcare Products. Any document in electronic form may be printed once. Copy and distribution of electronic documents and hardcopies is prohibited. Offenders will be liable for damages. All other rights are reserved.

Print No.: MR-000.816.17.05.02
Replaces: MR-000.816.17.04.02

English
Doc. Gen. Date: 07.11

Copyright

“© Siemens, 2004” refers to the copyright of a Siemens entity such as Siemens Aktiengesellschaft - Germany, Siemens Shenzhen Magnetic Resonance Ltd. - China, Siemens Shanghai Medical Equipment Ltd. - China, Siemens Medical Solutions USA Inc. - USA, Siemens Healthcare Diagnostics Inc. - USA and/or Siemens Healthcare Diagnostics Products GmbH - Germany.

Document Version

Siemens reserves the right to change its products and services at any time.

In addition, manuals are subject to change without notice. The hardcopy documents correspond to the version at the time of system delivery and/or printout. Versions to hardcopy documentation are not automatically distributed. Please contact your local Siemens office to order current version or refer to our website <http://www.healthcare.siemens.com>.

Disclaimer

Siemens provides this documentation “as is” without the assumption of any liability under any theory of law.

The installation and service of equipment described herein requires superior understanding of our equipment and may only be performed by qualified personnel who are specially trained for such installation and/or service.

1	General	7
	Introduction	7
	syngo MR system software	8
	Prerequisites	8
	MR 2004A Rev. B	8
	How to use this Register	8
	Replacing the HOST	8
	CCU - Firmware / Loadware / Testware	10
	Cameras	10
	Customer files	11
	Remote diagnostics	12
	Safety	13
	SAR	13
	Stimulation Monitor	13
	Access to the system	13
	Integrity of software back-up	13
	How to Navigate in this Manual	15
	Start points	15
	DICOM Configuration	16
2	Getting started and configuration pages	17
	General	17
	Passwords	17
	Prerequisites for the main console	17
	Prerequisites for the satellite console	17
	Site specific configuration	18
	Prerequisites MRC	18
	Prerequisites MSC	18
	Configuration pages for MRC	20
	System	20
	Database	24
	Service	24
	DICOM	28
	Import/Export	34
	EPR Interface	35
	External Devices	35
	Applications	35
	Measurement	38
	Configuration pages MRSC	40
	System	40
	Local host	41
	Database	42
	DICOM	45
	Import/Export	45
	External Devices	47
	Applications	47

3	Installing MRC / MRC+MRSC	51
	Overview of installing VA25A MRC/MRC+MRSC	51
	General	52
	Installing Software	53
	Prerequisites	53
	Installing MR 2004A Rev. B MRC / MRSC	54
	Keyboard settings and Digit grouping symbol for administrator	55
	Export MRC configuration data to MRSC	59
	Import configuration data from CD-R to MRSC	60
	Reboot	61
	Initial configuration after software installation MRC and MRSC	62
	Keyboard settings and Digit grouping symbol for meduser	62
	MRC	67
	Activate BIOS protection	67
	IFN (Image Filter New) Installation Work Steps	68
	Hints	68
	Installing IFN on MRC and MRSC	69
	Verifying correct IFN installation	70
	CV License (only valid for systems delivered before 01/2006)	70
	MRSC	71
	General	71
	MRSC Sharing Database with MRC	72
	Completed configuration	72
	Activate BIOS protection	72
	IFN (Image Filter New) Installation Work Steps	73
	Hints	73
	Installing IFN on MRC and MRSC	74
	Verifying correct IFN installation	75
	CV License (only valid for systems delivered before 01/2006)	75
4	Installing Software at the Imager	76
	Overview installing Imager SW VA25	76
	General	76
	Installing Software at the Imager "R630"	77
	Prerequisites	77
	Installation procedure	77
5	Site specific data	79
	Save site specific data	79
	Numaris4	79
	CustomerProtocols	79
	SW-Settings	79
	Prerequisites	79
	Settings that are not part of the backup/restore	80
	Backup (Save site specific data) on CD	80
	Backup (Save site specific data) on USB device	81
	Additional packages	82

	Restoring site specific data MR 2004A Rev. B on MR 2004A Rev. B.	84
	Prerequisites.	84
	Restore site-specific from CD.	84
	Restore site-specific from external USB device	86
	Settings that must be manually activated after restore	86
6	Configuration	88
	General.	88
	Prerequisites.	89
	Access Service Software	90
	Changes in the configuration	90
	Parallel Printer Installation	91
	Installing local parallel printer	91
	Enabling MRC to access printer share of MRSC	91
	Network Printer.	92
	Installing network printer.	92
	General Printer Properties	93
7	Remote connectivity	97
	General.	97
	Prerequisites.	98
	System Management	98
	Pro Active Service.	98
	Quality Assurance DB (QDB)	98
	Event DB (EDB)	98
	System Utilization	98
	Check / add route to SIEMENS Remote Services router	98
	Configure MAIL and FTP target (router systems).	100
	Testing connectivity from system to Rdiag server	102
	Send a test mail	102
	Test FTP login to the Remote Server	104
	AutoReport configuration	106
	Configuring QDB report	106
	Configuring EDB report.	107
	Configure System Utilization reporting	108
	System Utilization.	110
	System Management Installation	111
	Prerequisites.	111
	Routing	112
	Check connectivity to the System Management server	113
	Required ports for System Management (CA Unicenter) Communication	113
	Master.	114
	Troubleshooting and useful information	117
	System Management Installation failed	117
8	Appendix	120

Re installation of software fails	120
Formatting the system disk	120
System language and weight	122
Keyboard	123
Language Settings Japanese	124
Switch language from international (english-US) to japanese	124
Switch language from japanese to international (english-US)	127
Language Settings Chinese	130
Switch language from international (english-US) to chinese	130
Switch language from chinese to international (english-US)	133
Institution name	135
Regional settings	136
Date and Time syntax	136
System Manager	137
Restart Application	137
Copy SIEMENS Programs to the customer tree	138
IP Security	139
Adding host to IP security list	139
Deleting all IP security settings	139
Adding all configured hosts to IP security list	139
Context Vision	141
Requested data for cvlicense generation	141
Activate license for Context Vision	141
MR auto login	143
Access NT login	143
syngo startup / shutdown counter	144
How to terminate the syngo counter	144
Terminate deadlock of crashed Service UI	145
How to terminate the Service IU (Internet Explorer)	145
Expert - i	146
Measurement system configuration matrix	147
Dicom	153
Dicom and UNIX nodes	153
Dicom configuration data for SIEMENS products	153
9 Changes to previous version	156
10 List of Hazard IDs	157

Introduction

The installation and service of equipment described herein is to be performed by qualified personnel which is employed or authorized by Siemens Medical Solutions or one of its affiliates.

Assemblers and other personnel not employed by nor directly affiliated with nor authorized by Siemens Medical Solutions or one of its affiliates technical services are directed to contact the local office of Siemens Medical Solutions or its affiliates before attempting installation or service procedures.

syngo MR system software

Prerequisites

To configure the syngo-based MR system following training is required in combination with this Software manual:

- XP-administrator training
- Network training (TCP/IP)
- DICOM training
- syngo training / syngo MR training

MR 2004A Rev. B

The following software media must be delivered with the system: **CD for VA25A**

The MR 2004A Rev. B software includes the service software which is protected by a license mechanism.

NOTE

The "Release Notes syngo MR MR 2004A Rev. B" provides additional information regarding the MR 2004A Rev. B software.

MRC / MRSC R630

Operating system Windows XP Professional (syngo MR 2004A Rev. B)

Imager R630

Operating system Windows XP Professional (syngo MR 2004A Rev. B)

How to use this Register

Save time by reading this chapter first

This chapter contains reference, safety, and information regarding navigation in this manual.

Replacing the HOST

- The dongle from the defective host has to be installed at the new host due to the license mechanism.

- Applications licenses are generated in relation to the flexid of the 7, 9- (Dongle) with Census plus.
Exception: if a 8- Dongle is plugged, all application licenses are generated in relation to the MAC address of the first NIC (one exception is the Context Vision filter, it needs two licenses).
The low level license of CV filter is always generated in relation to the DONGLE.
The Service license is generated in relation to the DONGLE always (7, 9- and 8- dongles).
- IP settings from the own host might not be restored from the save set, if the old NIC PCB differs to the new NIC: Enter the local IP address manually in the service platform.
- If the Graphic card differs in between the original host and the new host: Perform the monitor calibration.

NOTE

All licenses for customer options are handled by TD HSC 1 (phone +49 9191-188080).

CCU - Firmware / Loadware / Testware

Cameras

NOTE

The use of non-released combinations of hardcopy cameras and BU products is prohibited and does not meet the requirements defined. For detailed information regarding the CAMERA configuration, refer to technical documentation "Camera Info Numaris 4 VA25A"

Customer files

User settings for Filming, Viewing, etc., are part of the save set package "**User-settings**" in SW-Settings_02.

The Customer Protocols are stored in the save set package "**CustomerProtocols**" (Customer protocols are part of the Numaris4 save set package as well).

Ask the customer to contact service for support after making major changes to user-defined protocols and user settings.

Remote diagnostics

Only IE 6.0 has been tested and released for remote service as service browser at syngo MR systems running with NUMARIS MR 2004A Rev. B.

The "SIEMENS Remote Services TM" will only support router connections for NEWLY delivered systems. See ([Siemens Remote Service, http://cs.med.siemens.de/cms/productinfo/default.asp?startfolder=e3371dfd0ace1e5181c7d36777ef3b4b](http://cs.med.siemens.de/cms/productinfo/default.asp?startfolder=e3371dfd0ace1e5181c7d36777ef3b4b)) for the most recent information regarding the router configuration..

(H WS CS - Customer Services -> For Service -> Product Information -> Siemens Remote Service)

The project manager is responsible to plan for a router connection to the Siemens Remote Diagnostic server. Regarding questions to your pre configured router from Siemens please contact mailto:srs_final@med.siemens.de or your remote diagnostic responsible.

Safety

SAR

The absorption of RF-energy by the patient is an integral part of every MR examination. The amount of energy absorbed has to be within the limits specified by the respective federal agency.

The Specific Absorption Rate (SAR) is different for each country. For this reason, you have to select the country code in the site-specific configuration corresponding to the local guidelines and regulations.

In addition, the temperature and relative humidity in the examination room affect the SAR values (only applies to countries subject to FDA and IEC guidelines, that is, all countries except Japan).

Stimulation Monitor

The magnetic field changes very quickly on Sonata and Quantum systems during gradient switching. If the electric field generated exceeds a specific threshold value, an electric current may be induced in the patient.

The Stimulation Monitor threshold is different for each country. For this reason, you have to select the country code in the site-specific configuration corresponding to the local guidelines and regulations.

The current induced may cause noticeable peripheral nerve stimulations and a feeling of discomfort for the patient.

Access to the system

The system software can only be accessed via password-protected accounts. This safety feature protects confidential (patient) data and reduces the possibility of (accidental) software manipulation by a third party.

The customer is provided with the passwords for the accounts required to operate and maintain the system.

Integrity of software back-up

The software has undergone extensive testing. However, this does not ensure that the software will run indefinitely without errors.

A backup (save set of site-specific data) should be created to avoid loss of system-specific data if the software has to be restored, e.g. after a system disk crash.

The save set will reduce system downtime to a minimum.

NOTE

Before you release the system to the customer, ensure that the image quality is within specifications and that all safety-relevant checks have been passed, e.g., the image orientation is correct.

Refer to Register "Tune-up" in the Technical Documents.

How to Navigate in this Manual

This manual includes information about various tasks. Please see the corresponding chapters listed in the table:

Start points

Getting started and configuration pages Your notes, workflow for initial configuration of new systems (software is pre-installed at the factory)	Refer to (Getting started and configuration pages / p. 17) 
Installing MRC / MRC+MRSC MRC or MRSC suffers from a defect that requires re-installation of the "syngo MR 2004A Rev. B"	Refer to (Installing MRC / MRC+MRSC / p. 51)  Refer to (Installing MRC / MRC+MRSC / p. 51) 
Installing Software at the Imager The Imager suffers from a defect that requires re-installation of the "syngo VA25A"	Refer to (Installing Software at the Imager / p. 76) 
Site specific data Save or Restore site-specific data, e.g. to save the configuration after making changes or restore site-specific data after software installation	Refer to (Site specific data / p. 79) 
Configuration Configuring / reconfigure a running system.	Refer to (Configuration / p. 88) 
Remote connectivity Configuring Quality Assurance Reporting and System Utilization reporting.	Refer to (Remote connectivity / p. 97) 
Appendix General information (e.g., BIOS settings and CV license, Institution name, Date, Time, Measurement system data).	

DICOM Configuration

The "CCC Team" provides a powerful configuration management tool see <http://cs.med.siemens.de/> "For Service", "DICOM", in the DICOM Homepage (CCC) navigate to "Interoperability DB", "Download", "Acquisition Tool".

General

This chapter describes the steps involved in configuring / reinstalling syngo-MR systems running VA25A.

MRC/MRSC are used as abbreviations for the Main/Satellite consoles, mrc/mrsc are used as host names.

This part of the table is always valid.

This part of the table is only valid for MRSC consoles
--

Passwords

NOTE

All passwords required to install and configure the system are available on the Med TD Intranet server at <http://www-td.med.siemens.de/>, "CS Knowledge Base", Speed Infos, search in "Speed Info" (except for the service key = variable password, which is handled by the districts).

- For more information regarding the service key (variable password), refer to the "Service Software" in TSG.

Prerequisites for the main console

- The VA25A software is installed at the MRC.

Prerequisites for the satellite console

- The VA25A software is installed at the MRC and MRSC.

NOTE

If all prerequisites are fulfilled, choose one of the suggested procedures and follow the instructions and the flowchart.

Site specific configuration

Recommendation: For paper versions, take out this page and use it for navigation.

Prerequisites MRC

- System software MR 2004A Rev. B is installed at the MRC

Prerequisites MSC

- Configuration is completed for the MRC

- System software MR 2004A Rev. B is installed at the MRSC

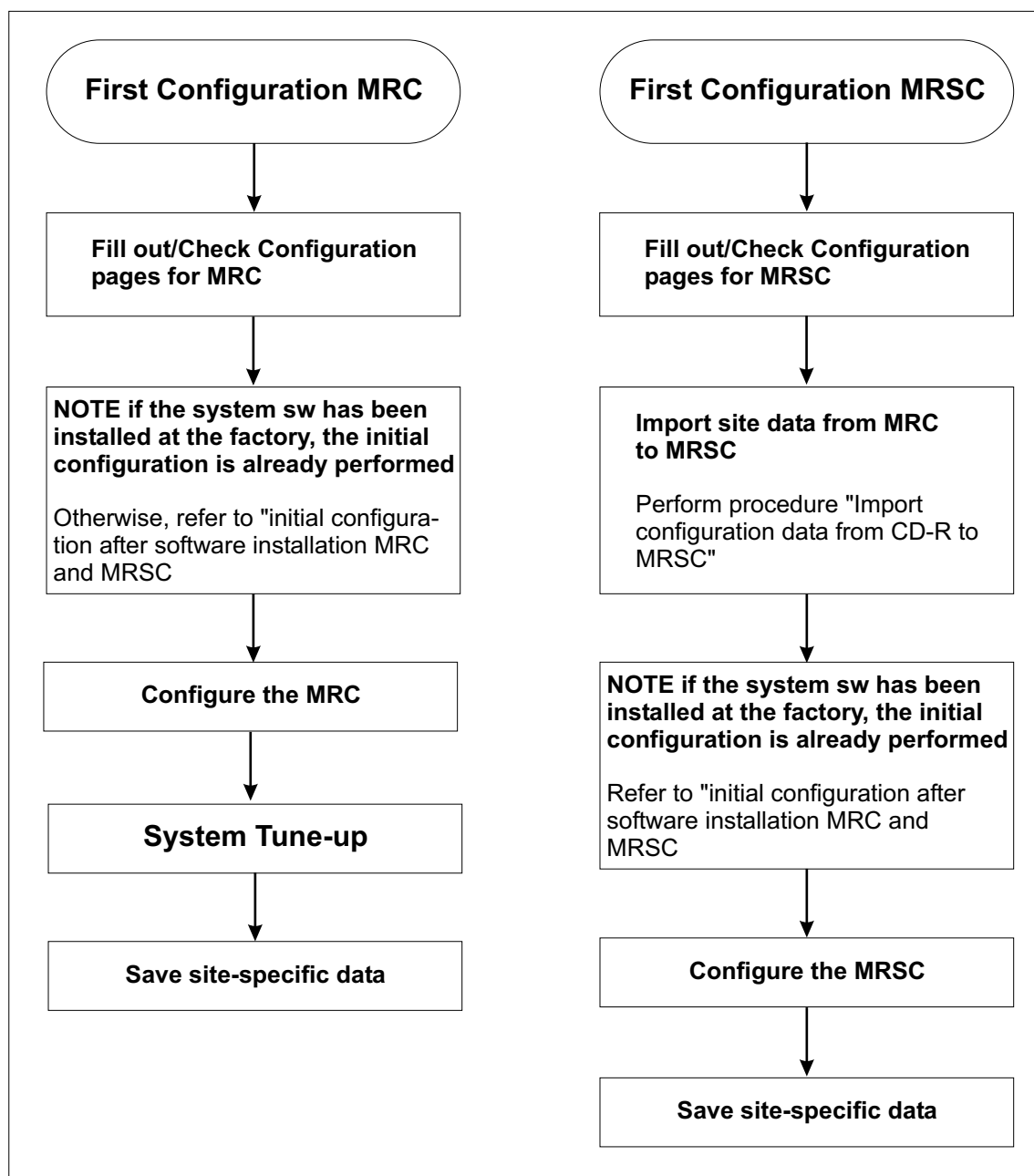


Fig. 1: First Configuration MRC/MRSC

Configuration pages for MRC

Before you start the configuration of the system, please fill out the following configuration pages. The order of the tables in the configuration pages corresponds to the order during “Full Configuration”.

Actual SW version	
Software version	Signature / Date

System

NOTE

If the system is connected to a hospital network, contact the network administrator regarding IP addresses and network settings.

list of system options	
Option / System type	Your notes
Attached to network	
Modem	
Main Console with Satellites or Satellite Console	
DICOM Print Devices	
DICOM Offline Devices	
DICOM Networking	
DICOM HIS/RIS	
Image Import/Export	
Paper Printer	
EPRI	

Local Host

MR System Specification		
Input field	Your notes	Note
Console Type / System Type		
System type e.g.: Harmony, Symphony, Sonata, Concerto, Har- monyExpert, SymphonyVi- sion, SonataVision, Allegra, or Trio		

Site Info		
Input field	Your notes	Note
System		
Serial No.		
Maintenance Unit	mr	
Department		
Service		
Site Identification No		
Service Center Phone		
Customer and Address		
Name		
Id		
Hospital		
Street		
Street No		
Zip Code		
Phone No		
City		
District		
Country		
Time zone		
		Pull-down menu

Site Info		
Input field	Your notes	Note
AutoDaylightTime		
(Day light saving)		Check box

NOTE

The host name must be unique in the network!

TCP/IP LAN		
Input field	Your notes	Note
Identification		
Computer name	Recommended syntax: e.g. MRC22000 or MRSC22000	(MRC and system serial number)
		Host name for MRC Use only letters from a to z (numbers are allowed as well, but they may not be used for the first character of the host name)
Adapter (Network adapter)		
select name		Pull-down menu lets you select from a predefined list
IP Address		
IP address (local host)		Delete the defaults in the configuration menu for the first Network PCB (only one entry is supported in this version)
Subnet mask		
gateways		
NOTE: Use default addresses for the second Network PCB 192.168.2.1 (NIC for measurement system network MPCU and Imager)		
WINS		
Primary WINS server		
Secondary WINS server		
Enable DNS for Windows Resolution		Select box
Enable LMHOSTS Lookup		Select box

TCP/IP LAN		
Input field	Your notes	Note
DNS		
Domain		
1 DNS service search order		Several DNS nodes can be defined
2 DNS service search order		
3 DNS service search order		

Users		
Input field	Your notes	Note
Account for UI	Leave default Default is meduser	Changing the password / account will result in problems during manual adjustments
Default Password / or give	Leave default password	
Account for PPP (Account name is: PPP)		
Default Password / or give		Account name: PPP This account is required for remote service sessions via modem. You must agree on a unique PPP password for your system together with the responsible USC/HSC Rdiag server administrator
Enable Autologin		Autostart of the MR application

Monitor Type settings		
Input field	Your notes	Note
Monitor Type		
color / greyscale		For color LCD monitors select: Appearance Siemens (True Color)

Database

Satellite Console (host name and IP-address of the MRSC)		
Input field	Your notes	Note
Host Properties (MRSC)		
Select Host		Enter the IP address of the MRSC and host name Only LAN is supported!

Service

Autoreport (mail target settings for service autoreporting tool)		
Input field	Your notes	Note
setting of LAN connection of remote mail server		
Select Host		Several target hosts can be set Note the logical name
Host Properties (remote server host)		
1 Host Name		Note the host names and IP addresses
1 TCP/IP address		
2 Host Name		
2 TCP/IP address		
3 Host Name		
3 TCP/IP address		
4 Host Name		
4 TCP/IP address		

Mail		
Input field	Your notes	Note
Host Properties		
Select Host		Enter the IP address of the Rdiag server

Mail		
Input field	Your notes	Note
1 Host Name		Note the host names and IP addresses
1 TCP/IP address		
2 Host Name		
2 TCP/IP address		
Mail Server Host Selection		
email address		Enter the E-mail address of your local USC
email address in Remote Service Access Control window		Check-box (If this check-box is set the customer function Mail in the "Remote Access Control window" will be sent to the above e-mail address)

FTP (ftp target settings for service file&image transfer tool)		
Input field	Your notes	Note
Select Host		Several target hosts can be set
		Note the logical name
Host Properties (Remote server host)		
1 Host Name		Note the host names and IP addresses
1 TCP/IP address		
2 Host Name		
2 TCP/IP address		
Logical Name		
FTP Host Selection		Pull-down menu
FTP User		
FTP Password		(Confirm password at the second row)

AutoTransfers		
Input field	Your notes	Note
activate/deactivate automatic Transfers (initiated by the autoreport tool)		
activate automatic Transfers		Must be enabled that Auto Reports can be sent automatically by the system. system utilization will only work if check box is enabled.

Eventlog		
Input field	Your notes	Note
Configuration of display parameters for the Eventlog		
Log only Service messages		Check box
Display Service messages		Check box
Display sequence number		Check box

Backup/Restore		
Input field	Your notes	Note
package		
Additional package 1		The following packages are pre selectable by default: Numaris4 SW-Settings CustomerProtocols Diagnostic ExamDB Additional packages might be included for re
Additional package 2		
Additional package 3		
Additional package 4		
Additional package 5		
Additional package 6		
groups		

Backup/Restore		
Input field	Your notes	Note
files & directories		
file & dir		
backup pre		
backup post		
restore pre		
restore post		

Licensing		
Input field	Your notes	Note
Input/Change of Licenses		
License Key		Keep the list with licences separated from this documentation
Hardware ID Type	Dongle	
Hardware ID		Unique dongle number (ID)

Fonts		
Input field	Your notes	Note
Input/Change of Licenses		

Fonts		
Input field	Your notes	Note

DICOM

The "CCC Team" provides a powerful configuration management tool (see "DICOM", CCC Connecting Systems, "Databases", "Interoperability db", "Download", "Acquisition Tool")

(Acquisition Tool,
https://intranet.healthcare.siemens.com/cms/IKM/en/Med-wideTopics/Connectivity_Competence_Center/Databases/Interoperability_DB/Pages/Acquisition_Tool.aspx)

NOTE

Ask your hospital network administrator for the correct DICOM nodes.

General		
Input field	Your notes	Note
Local application entity titles		
HIS / RIS		(SCU)
Study Transfer		(SCP + SCU)
Print		(SCU)
Simultaneous DICOM associations		
maximum PDU size		(Packet Data Unit size)
time-out values	use default settings	
implicit raw data transfer		Check-box

character set (configuration of the DICOM default character set used when registering new patients)

Input field	Your notes	Note
DICOM default character set		

character set (configuration of the DICOM default character set used when registering new patients)		
Input field	Your notes	Note
	Set radio button to "allow multiple character sets"	ISO 2022 IR100 Japan Only: To support Japanese character set, select ISO 2022 IR13
Select character set		Pull-down menu

Offline devices (Setting of DICOM archive devices on local host)		
Input field	Your notes	Note
Select Logical Name		
Logical name		Several targets can be set e.g CDRW MOD etc.
Logical name		
Logical name		For MOD drives: set the SPI file system flag
Logical name		

Network Nodes setting of LAN connection of remote DICOM transfer target)		
Input field	Your notes	Note
Host Properties		
Select Host	1	Several target hosts can be set Note the logical name
	2	
	3	
	4	
	5	
	6	
Host Properties (remote server host)		

Network Nodes setting of LAN connection of remote DICOM transfer target)		
Input field	Your notes	Note
Host Name	TCP/IP address	Note the host names and IP addresses
1		
2		
3		
4		
5		
6		

Network Nodes (definition of application entity titles and ports of DICOM storage SCPs)		
		Note
General Node Properties		
Logical Name	(Default node / syngo based)	(syngo based yes/no, first default node or second default node etc.)
1.	1.	
2.	2.	
3.	3.	
4.	4.	
5.	5.	
6.	6.	
Application Entity Properties		
AE Title	Port Number	
1.	1.	
2.	2.	
3.	3.	
4.	4.	
5.	5.	
6.	6.	
Supported DICOM services		
Storage, Storage, Commitment, Query, and/or Retrieve		

Network Nodes (definition of application entity titles and ports of DICOM storage SCPs)		
		Note
1.		If you are unsure which services are needed, please select all
2.		
3.		
4.		
5.		
6.		

Print Devices (definition of application entity titles and ports of DICOM storage SCPs)		
		Note
General Node Properties		
Logical Name	(Default node / syngo based)	(syngo based yes/no, first default node or second default node etc.)
1.	1.	
2.	2.	
3.	3.	
4.	4.	
5.	5.	
6.	6.	
Application Entity Properties		
AE Title	Port Number	
1.	1.	
2.	2.	
3.	3.	
4.	4.	
5.	5.	
6.	6.	
Supported DICOM services		
Print		

hc overview (overview of hardcopy devices)		
Input field	Your notes	Note
Select logical name 1		
Name		Remember to mark the Default Camera
Type		
Comment		
Location		
Hold printed filmjobs		(Pulldown menu)
Network access		
DICOM access		
Select logical name 2		
Name		
Type		
Comment		
Location		
Hold printed filmjobs		(Pulldown menu)
Network access		
DICOM access		
Select logical name 3		
Name		
Type		
Comment		
Location		
Hold printed filmjobs		(Pulldown menu)
Network access		
DICOM access		

NOTE

The import of LUT files is only required in certain cases.

Image type / Max Density		
Input field	Your notes	Note
LUT	MR1	Use defaults (NEW with MR 2004A Rev. B LUT with 16 bit are available Presentation LUT is supported.
Density	3.00	

HIS/RIS nodes (setting of LAN connection to remote HIS/RIS station)		
Input field	Your notes	Note
Select Host		Several target hosts can be set Note the logical name
Host Properties (remote server host)		
1 Host Name		Note the host names and IP addresses
1 TCP/IP address		
2 Host Name		
2 TCP/IP address		
3 Host Name		
3 TCP/IP address		
4 Host Name		
4 TCP/IP address		
Select Logical Name		
Logical name		Several targets can be set e.g. CDRW MOD etc. For MOD drives: set the SPI file system flag
Logical name		
Logical name		
Logical name		
definition of application entity title and port of DICOM report or worklist SCP		
General Node Properties		
Logical Name		
Host		

HIS/RIS nodes (setting of LAN connection to remote HIS/RIS station)		
Input field	Your notes	Note
Application Entity Properties		
AE Title		
Port Number		
Supported DICOM services		
Basic Worklist		Check box
Performed Procedure Step		Check box

Import/Export

directories (setting of LAN connection to remote computer)		
Input field	Your notes	Note
Host Properties		
Select Host	1	Several target hosts can be set Note the logical name
	2	
	3	
Host Properties (remote server host)		
Host Name	TCP/IP address	Note the host names and IP addresses
1		
2		
3		
Select logical HD device		
Logical Name	Path	Drive
1		
2		
3		
Remote system		

directories (setting of LAN connection to remote computer)		
Input field	Your notes	Note
Host	Share	
1		
2		
3		

EPR Interface

Server		
Input field	Your notes	Note
Location of EPRServer		http://123.456.789.10/webEPR
Type of EPRServer	SHS	

External Devices

Paper Printer		
Input field / Your notes		
Patient		
Letter small:	A3:	A4:
A3:	B4:	B5:
A5:	Legal:	

Applications

Pat. Registration		
Input field		Your notes
Displayed Patient Registration Dialog Elements		
Patient	Study	
Name Fields	Study List	
Age	Accession No	
Height	Requested Procedure	

Pat. Registration		
Input field		Your notes
Weight	Request Id	
Address	Study Comment	
Other Name(s)		
Other ID(s)		
Ethnic Group		
Military Rank		
Details	Visit	
Additional Info	Referring Physician	
Institution	Requesting Physician	
	Admitting Diagnosis	
	Ward	
	Admission Id	
Patient Registration Dialog Invocation		
implicitly in scheduler		

Viewer		
Input field		Your notes
Viewer configuration settings		
Close all studies when new one is loaded		
Close all series when new one is loaded		

CorRea		
Input field		Your notes
Displayed CorRea Dialog Elements		

CorRea		
Input field		Your notes
Patient	Instance	
Name Fields	Number	
Other Name(s)	Date	
Other ID(s)	Time	
Ethnic Group	Contrast Medium	
Military Rank	Comment	
Height	Calibration Image	
Weight	Image Laterality	
Address	Series	
Additional Info	Request ID	
Hospital	Modality	
Referring Physician	Patient Position	
Requesting Physician	Requesting Physician	
Admitting Diagnose	Institution Address	
Ward	Performing Physician	
Admission ID	Operator	
Study	Body Part Examined	
Accession Number	Number	
ID	Description	
Description	Product Name	
Comments	Date	
Date	Time	
Time	Laterality	

MPPS		
Input field		Your notes
Displayed MPPS Dialog Elements		
Study	Study ID	
Study Date	Study Time	
Accession No.	Description	

MPPS		
Input field		Your notes
Modality	Station	
Status	Type	
End Date	End Time	
Dose Dialog	Comments	

Worklist Results		
Input field		Your notes
Displayed HISRIS worklist result attributes		
Date of birth	Sex	
Accession Number	Admission ID	
Referring Physician	Performing Physician	
Requesting Physician	Ward	

Measurement

System type		
Input field	Your notes	Note
Patient		
MPCU type:		MPCU/MPCU-2/MPCU-3:
# of RX4 boards		1 or 2
# of RF receiver		1, 2 or 4
gradient amplifier		Type of gradient AMP
gradient coil		
Stimulation Monitor		
SAR guideline		
Advanced SHIM option		
PMU option		
Unit type		(stationary / mobile)
MRIR settings		
# of Harddisks		1 or 4

Only valid for MAGNETOM Concerto (007)		
Input field	Your notes	Note
Easy site	Yes / No	Sites without external interferences, e.g. railway etc. Refer to site measurement protocol
EFI	none / installed	
mains frequency (EFI)	50 Hz / 60 Hz	Mains frequency on site
railway mains frequency (EFI)	16,7 Hz or 25 Hz	Only valid if a railway is close to the site

Configuration pages MRSC

Before you start the configuration of the system, please fill out the following configuration pages. The tables in the configuration pages are listed in the same order as they are requested during “Full Configuration” in install.

NOTE

For the MRSC, only the configuration of the following items is important:

- Computer name
- IP address (local host)

Actual SW version	
Software version	Signature / Date

System

list of system options	
Option / System type	Your notes
Attached to Network	MUST be checked
Main Console with Satellites or Satellite Console	MUST be checked
DICOM Print Devices	
DICOM Offline Devices	
DICOM Networking	
DICOM HIS/RIS	
Paper Printer	
Paper Printer	
EPRI	

Local host

MR System Specification		
Input field	Your notes	Note
Console Type		
Satellite Console	Satellite Console	

NOTE

If the system is connected to a hospital network, contact the network administrator regarding IP addresses and network settings.

NOTE

The host name must be unique in the network!

TCP/IP LAN		
Input field	Your notes	Note
Identification		
	Recommended syntax: e.g. MRSC22000	(MRSC and system serial number)
Computer name		Host name for MRSC Use only letters from a to z (numbers are allowed as well, but they may not be used for the first character of the host name)
Adapter (network adapter)		
select name		Pull-down menu lets you select from a predefined list
IP Address		
IP address (local host)		Delete defaults in the configuration menu for the Network PCB (only one entry is supported in this version)
Subnet mask		
gateways		
WINS		
Primary WINS server		
Secondary WINS server		

TCP/IP LAN		
Input field	Your notes	Note
Enable DNS for Windows Resolution		Select box
Enable LMHOSTS Lookup		Select box
DNS		
Domain		
1 DNS service search order		Several DNS nodes are definable
2 DNS service search order		
3 DNS service search order		

Users (Configuration of user accounts which will be needed for operation)		
Input field	Your notes	Note
Account for UI		Default is meduser
Default Password / or give	Leave the default Default is meduser	Changing the password / account will result in problems during manual adjustments
Enable Autologin	Leave default password	

Monitor Type settings		
Input field	Your notes	Note
Monitor Type		
color / greyscale		For color LCD monitors select: Appearance Siemens (True Color)

Database

Satellite Console		
Input field	Your notes	Note
Host Properties		Enter the hostname and IP address of the MRC (Main console) Only LAN is supported
Host Name		
TCP/IP address		
Console Type		

Satellite Console		
Input field	Your notes	Note
	Satellite Console	
Database Properties		
Name of Local Main DB		
Name of Local Scheduler DB		
Database Offset of local DBs [100..65500]	101	For multiple satellite consoles please set the offset values with proper spacing e.g. 100 for the first console and 200 for the second one. The main data base keeps offset 0

Eventlog		
Input field	Your notes	Note
Configuration of display parameters for the Eventlog		
Log only Service messages		Check box
Display Service messages		Check box
Display sequence number		Check box

Licensing		
Input field	Your notes	Note
Input/Change of Licenses		
License Key		Keep the list of licenses separate from this docu- mentation
Hardware ID Type	Dongle	
Hardware ID		Unique dongle number (ID)

Fonts		
Input field	Your notes	Note
Input/Change of Licenses		

Fonts		
Input field	Your notes	Note

NOTE

Following part of the configuration for the MRSC can be comfortably imported with the procedure [\(Export MRC configuration data to MRSC / p. 59\)](#).

Following configuration items

- LocalHost\Users
- TCPIP\LAN-Settings
- Database\SateliteConsole
- Database\DBOffset
- Applications\WorklistResults
- Applications\MPPS
- Applications\CorRea
- Applications\Viewer
- Applications\PatientRegistration
- EPRI\Server
- DICOM\ImportExportDirectories
- DICOM\HISRISNodes
- DICOM\LUTFiles
- DICOM\PrintDevices
- DICOM\NetworkNodes
- DICOM\OfflineDevices
- DICOM\General
- ExternalDevices\PaperPrinter
- Service\Licensing
- Service\Eventlog
- Service\AutoTransfers
- Service\FTP
- Service\Mail
- LocalHost\MonitorType
- LocalHost\SiteInfo

DICOM

Refer to your notes for MRC

Import/Export

NOTE

Use same settings as for MRC (Only in some cases the customer needs different setting in-between MRC and MRSC)

TCP/IP LAN		
Input field	Your notes	Note
Identification		
	Recommended syntax: e.g. MRSC22000	(MRSC and system serial number)
Computer name		Host name for MRSC Use only letters from a to z (numbers are allowed as well, but they may not be used for the first character of the host name)
Adapter (network adapter)		
select name		Pull-down menu lets you select from a predefined list
IP Address		
IP address (local host)		Delete defaults in the configuration menu for the Network PCB (only one entry is supported in this version)
Subnet mask		
gateways		
WINS		
Primary WINS server		
Secondary WINS server		
Enable DNS for Windows Resolution		Select box
Enable LMHOSTS Lookup		Select box
DNS		
Domain		

TCP/IP LAN		
Input field	Your notes	Note
1 DNS service search order		Several DNS nodes are definable
2 DNS service search order		
3 DNS service search order		

directories (setting of LAN connection to remote computer)		
Input field	Your notes	Note
Host Properties		
Select Host	1	Several target hosts can be set Note the logical name
	2	
	3	
Host Properties (remote server host)		
Host Name	TCP/IP address	Note the host names and IP addresses
1		
2		
3		
Select logical HD device		
Logical Name	Path	Drive
1		
2		
3		
Remote system		
Host	Share	
1		
2		
3		

External Devices

Paper Printer		
Input field / Your notes		
Patient		
Letter small:	A3:	A4:
A3:	B4:	B5:
A5:	Legal:	

Applications

NOTE

Use same settings as for MRC (Only in some cases the customer needs different setting in-between MRC and MRSC)

Pat. Registration		
Input field		Your notes
Displayed Patient Registration Dialog Elements		
Patient	Study	
Name Fields	Study List	
Age	Accession No	
Height	Requested Procedure	
Weight	Request Id	
Address	Study Comment	
Other Name(s)		
Other ID(s)		
Ethnic Group		
Military Rank		
Details	Visit	
Additional Info	Referring Physician	
Institution	Requesting Physician	
	Admitting Diagnosis	
	Ward	
	Admission Id	

Pat. Registration	
Input field	Your notes
Patient Registration Dialog Invocation	
implicitly in scheduler	

Viewer	
Input field	Your notes
Viewer configuration settings	
Close all studies when new one is loaded	
Close all series when new one is loaded	

CorRea	
Input field	Your notes
Displayed CorRea Dialog Elements	

CorRea		
Input field		Your notes
Patient	Instance	
Name Fields	Number	
Other Name(s)	Date	
Other ID(s)	Time	
Ethnic Group	Contrast Medium	
Military Rank	Comment	
Height	Calibration Image	
Weight	Image Laterality	
Address	Series	
Additional Info	Request ID	
Hospital	Modality	
Referring Physician	Patient Position	
Requesting Physician	Requesting Physician	
Admitting Diagnose	Institution Address	
Ward	Performing Physician	
Admission ID	Operator	
Study	Body Part Examined	
Accession Number	Number	
ID	Description	
Description	Product Name	
Comments	Date	
Date	Time	
Time	Laterality	

MPPS		
Input field		Your notes
Displayed MPPS Dialog Elements		
Study	Study ID	
Study Date	Study Time	
Accession No.	Description	

MPPS		
Input field		Your notes
Modality	Station	
Status	Type	
End Date	End Time	
Dose Dialog	Comments	

Worklist Results		
Input field		Your notes
Displayed HISRIS worklist result attributes		
Date of birth	Sex	
Accession Number	Admission ID	
Referring Physician	Performing Physician	
Requesting Physician	Ward	

Overview of installing VA25A MRC/MRC+MRSC

Recommendation: For paper versions, take out this page and use it for navigation.

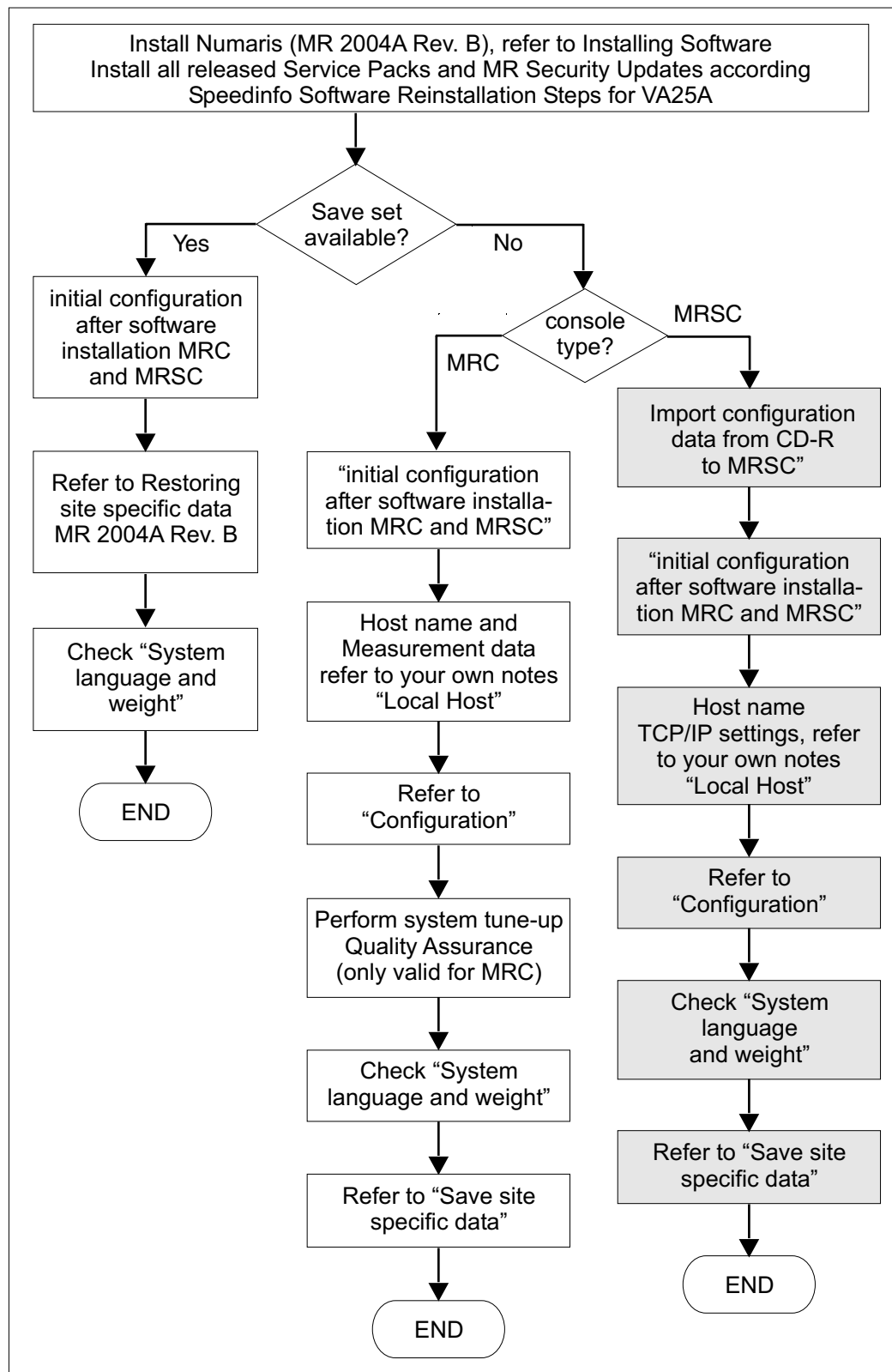


Fig. 2: Overview of installing VA25A MRC/MRC + MRSC

General

The following chapter describes the software installation of VA25A on an empty system disk at the MRC and MRSC for MAGNETOM systems.

NOTE

Make sure that you have all customer licenses as well as backups of site-specific data before reinstalling the software.

CV license filter - see ([Requested data for cvlicense generation / p. 141](#)).

NOTE

The Shim files are part of the dynamic data.

NOTE

In case of host replacement the local TCP IP settings are not restored.

NOTE

If the system has suffered from HD crash refer to ([MR-000.840.11 / Software](#)) TSG for troubleshooting information before you continue with the SW installation.

Table entries highlighted in gray are valid only for the MRSC installation:

This part of the table is always valid.

This part of the table is only valid for MRSC consoles

Installing Software

Prerequisites

- BIOS CD "MR059 BIOS CD VA25E with part no.: 100 49 286P9901"
- CD: 1/6 (International OS) 2/6 (XP Professional, syngo VD20N, CD 3/6 (NUMARIS/4 VA25A) is available.
- **(In case of MRC / MRSC host replacement including dongle. The "Dallas dongle" (Black) the application licenses are generated in relation to the Ethernet ID of the NIC. New application licenses must be generated in relation to the new Host)**

NOTE

The application licenses are part of the save set (back-up)
Only if the back-up is missing you can use the original delivered Licence keys on floppy or CD disk.

NOTE

The Service Licence is required for the configuration.

- Disconnect the LAN cable(s)

NOTE

Consoles with the optional Pioneer 1.7 MOD, the power must be switched on for the MOD drive during the software installation.

- Enter the BIOS

R630
press F2
Enter password
Press F9 -Y
Insert BIOS CD "MR059 BIOS CD VA25E withpart no.: 100 49 286 P9901"
Boot Device Priority
move CD-ROM to the top (with the + key)
Press F10 - Y

- Wait until: "System Ready for Software Installation" is displayed.
- Remove the BIOS CD

Installing MR 2004A Rev. B MRC / MRSC

- Switch on the console to be installed (MRC or MRSC).
 - Insert the **(MR 2004A 1/6 International OS or CD5/6: Chinese OS or CD6/6: Japanese OS)** into the CD-ROM drive. **If the SW re installation fails e.g. due to a interrupt during the initial installation part of CD 1/6 refer to [\(Re installation of software fails / p. 120\)](#)**

Output	Input	Note
SELECTION WINDOW continue is highlighted	↵	Press any key to continue
	Automatic reboot ...	
SELECTION MENU Distribution type: Customer ... Create Partitions on ... CONTINUE INSTALLATION		Distribution Select CONTINUE INSTALLATION with the arrow key Press ENTER (↵)
Confirm that "All DATA on selected harddisk) will be lost"	↵	Press ENTER (Installation time ~7 minutes)
INFORMATION WINDOW Please insert the next installation media labeled: syngo MR 2004A CD 2/6 Press any key to continue .. - Copying in progress		Please insert CD 2/6 Press ENTER (↵) (Installation time ~4 minutes)
INFORMATION WINDOW Please insert the next installation media labeled: syngo MR 2004A CD 3/6 Press any key to continue .. - Copying in progress		Please insert CD 3/6 Press ENTER (↵) (Installation time ~4 minutes)
The CD-ROM drive opens the carriage automatically for CD 3/6		Remove the CD from the CD-Rom drive

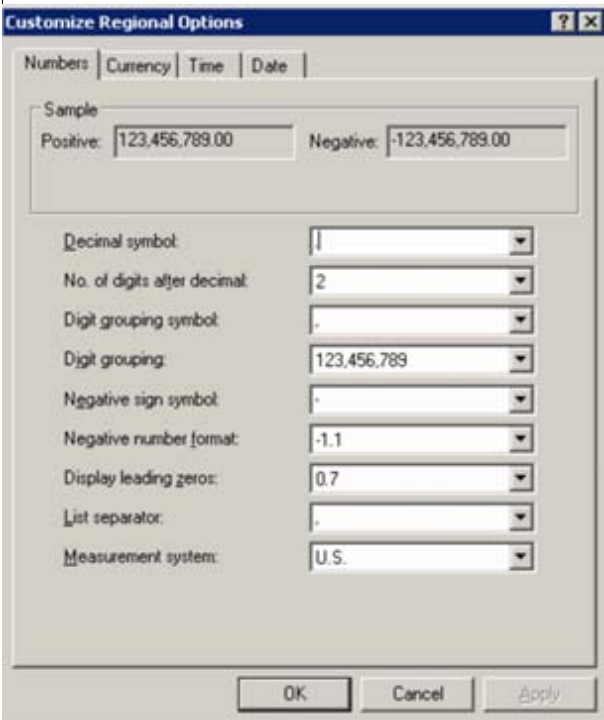
Output	Input	Note
Automatical reboot Windows XP Professional Setup Please wait while Setup copies files to your hard disk	Do not interfere the automatic installation / configuration. Windows XP Professional Setup will follow Several re-boots will follow Estimated time for the automatic configuration ~72 minutes	
After installation is complete, the syngo MR 2004A Logon window is displayed.		
Install all released Service Packs and MR Security Updates according Speedinfo Software Reinstallation Steps for VA25A		

Keyboard settings and Digit grouping symbol for administrator

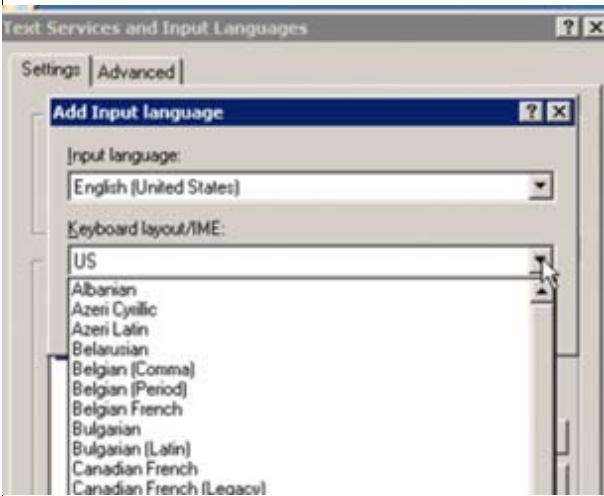
Following procedure is only to be performed for systems with **non-US** keyboard.

- Login as user **administrator**

Output	Input	Note
	Click START -> Settings -> Control Panel	Setting the local keyboard type
 Date, Time, Language, and Regional Options	Click Date, time, Language, and Regional options at the bottom right in the control panel	Only to be performed for systems with non-US keyboard
 Regional and Language Options	Click Regional and Language Options at the bottom right in the control panel	

Output	Input	Note
 The 'Regional and Language Options' dialog box is shown with the 'Regional Options' tab selected. Under 'Standards and formats', the 'Customize' button is highlighted with a mouse cursor. Below this, sample values for Number, Currency, Time, Short date, and Long date are displayed. The 'Location' section shows 'United States' selected.	Click Customize	
<i>Fig. 3:</i>  The 'Customize Regional Options' dialog box is shown with the 'Numbers' tab selected. The 'Sample' section shows 'Positive: 123,456,789.00' and 'Negative: -123,456,789.00'. Below this, various number formatting options are listed with dropdown menus: Decimal symbol (.), No. of digits after decimal (2), Digit grouping symbol (,), Digit grouping (123,456,789), Negative sign symbol (-), Negative number format (-1.1), Display leading zeros (0.7), List separator (.), and Measurement system (U.S.).	Select the Numbers register Set "Decimal Symbol" to "." (= Full stop) and "Digit grouping symbol" to "," (= Comma). Click OK	
<i>Fig. 4: Decimal / Digit grouping symbols</i>		

Output	Input	Note
 Fig. 5:	Click at the Languages, register; Click the Details Button Click Add	

Output	Input	Note
 Fig. 6:	Click pull-down menu "Keyboard layout/IME" Select input language related to your keyboard layout Click OK	
 Fig. 7:	Click pull-down menu "Select one of the installed input languages" Click Apply Click OK	The key-board settings will be valid only for user administrator

NOTE

If the system consists of MRC and MRSC check that both consoles are allowed to access each other through the ipsec filter.

The syngo-MR systems are protected with a service "ipsec".

Other network nodes, like MRSC or MRC and DICOM nodes, can only access the DICOM services / DB access that your console provides if they are included in the "Generic Transport Filter" of ipsec.

To list nodes that are allowed to access your console, click **Reports**; select **IP Security** below **SiteInfo**. To add a node manually into the IP security list, use following procedure:

At the "MRC", simultaneously press **Ctrl - Esc**. Select **Start, Advanced User**; enter the advanced user password; click **OK**

Simultaneously press **Ctrl - Esc**, click **Start, Run cmd** enter the command: **runas /user:administrator cmd** enter the password

In the newly started cmd tool you have administrator rights.

(replace the address "146.254.222.333" in the following **example** with the IP address of your remote node)

Input in the cmd window	Output
ipseccmd -f (0+146.254.222.333)	The command completed successfully

What to do next?		
Your Console type?		
MRC	MRSC with a valid save set	MRSC WITHOUT save set
Re-boot the console see (Reboot / p. 61) 		Keep login as administrator see (Export MRC configuration data to MRSC / p. 59) 

Export MRC configuration data to MRSC

NOTE

If you have a valid save set for the MRSC skip this procedure. Continue with ([Reboot / p. 61](#))

It makes sense only to perform the save at CD-R when you have completed the configuration of Customer UI: syngo, DICOM: Hardcopy devices and Network nodes at the MRC.

- Login as **administrator** at the **MRC**
- Insert an empty CD-R disk into the CD-R drive of the **MRC**

Output	Input	Note
	Click START Select RUN	Click START in the lower left corner Select RUN
	cmd Click OK	
	mrscsave.bat ↵ Press Enter	

Import configuration data from CD-R to MRSC

Only syngo UI settings, DICOM configuration data and TCP/IP addresses from remote nodes are copied with this procedure.

Systems consisting of Main and Satellite consoles

- You are logged in / login as **administrator** at the **MRSC**
- Insert CD-R disk into the CD-ROM drive of the **MRSC**, containing partial configuration data from the MRC (saved with the procedure "Export MRC configuration data to MRSC")

Output	Input	Note
	Click START Select RUN	Click START in the lower left corner Select RUN
	cmd Click OK	
	mrsc_default.bat ↵	

NOTE

The MRSC will use the database of the MRC for export/import to remote DICOM nodes for shared DB configuration.

Reboot

Output	Input	Note
	Close all open windows	
 Fig. 8:	Click Start -> Shut Down Select Restart Click OK	Reboot the host

What to do next?

(Initial configuration after software installation MRC and MRSC / p. 62)



Initial configuration after software installation MRC and MRSC

Leave the customer LAN disconnected

NOTE

The "Ping" command does not work during the initial configuration since the LAN is disconnected.

- Login as user **meduser**

NOTE

Before you enter into the displayed service platform, check/set keyboard settings.

Keyboard settings and Digit grouping symbol for meduser

Following procedure is only to be performed for systems with **non-US** keyboard.

- Login as user **meduser**

Output	Input	Note
	Click START -> Settings -> Control Panel	Setting the local keyboard type
 Date, Time, Language, and Regional Options	Click Date, time, Language, and Regional options at the bottom right in the control panel	Only to be performed for systems with non-US keyboard
 Regional and Language Options	Click Regional and Language Options at the bottom right in the control panel	

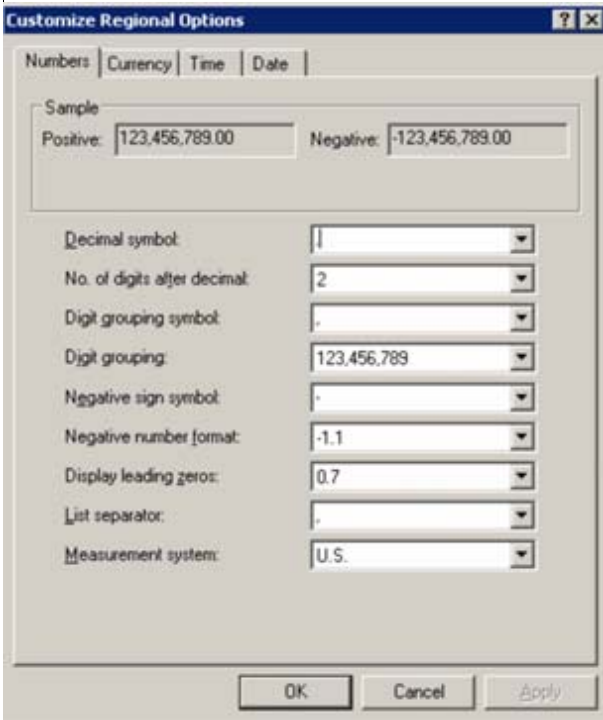
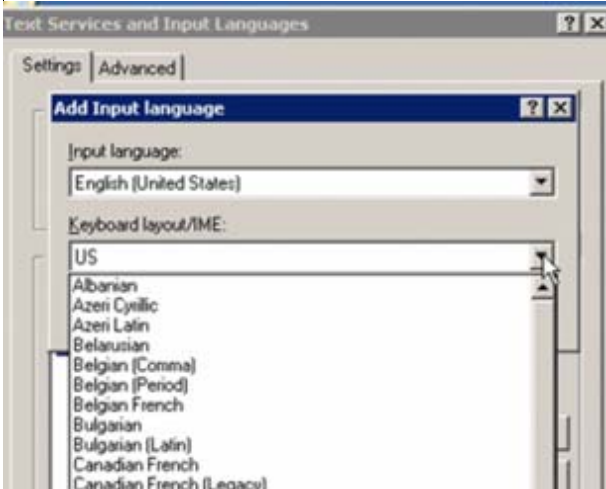
Output	Input	Note
 The 'Regional and Language Options' dialog box is shown with the 'Regional Options' tab selected. Under 'Standards and formats', the 'Customize' button is highlighted. The 'Location' is set to 'United States'.	Click Customize	
 The 'Customize Regional Options' dialog box is shown with the 'Numbers' tab selected. The 'Decimal symbol' is set to '.' and the 'Digit grouping symbol' is set to ','. The 'List separator' is also set to ','.	Select the Numbers register Set "Decimal Symbol" to "." (= Full stop) and "Digit grouping symbol" to "," (= Comma). Click OK	

Fig. 9:

Fig. 10: Decimal / Digit grouping symbols

Output	Input	Note
 Fig. 11:	Click at the Registrar Languages, Click at Button Details Click Add	

Output	Input	Note
 Fig. 12:	Click pull-down menu "Keyboard layout/IME" Select input language related to your keyboard layout Click OK	
 Fig. 13:	Click pull-down menu "Select one of the installed input languages" Click Apply Click OK	The keyboard settings will be valid only for meduser

NOTE

Ignore message "Could not retrieve security data"

NOTE

Keep login during the complete initial configuration, **DO NOT** re-boot!

Output	Input
	Click the Local Service link to start configuration
 The image shows the Siemens Service Software Authentication screen. At the top is the Siemens logo. Below it, the text 'Authentication' is in a black bar. The main title is 'Service Software'. Below that, it says 'Erlangen', 'Harmony MR 11002', and 'NUMARIS/4'. A prompt 'Please enter password' is followed by a text box containing '9CABDAF5B76323' and an 'OK' button. At the bottom, there is a checkbox labeled 'Set as Default'. <i>Fig. 14:</i>	Remember to click Set as default Enter the Service key (Level 7 key is required) Click OK

What to do next?	
Site specific data available? YES	Site specific data available? NO
See "Restoring site specific data MR 2004A Rev. B on MR 2004A Rev. B" on Page 5 - 4  The image shows two buttons: 'Backup' in blue and 'Restore' in red.	Perform full configuration (Remember to import the customer licences) Refer to your own NOTES Chapter 2 "Getting started and configuration pages" (It's recommended to use the configuration pages for navigation purpose where you are in the configuration procedure) NOTE Information regarding known problems during the Initial configuration are described here separately for MRC and MRSC.  The image shows a button labeled 'Config..' in blue.
	Console type MRC: See (MRC / p. 67) Console type MRSC: See (MRSC / p. 71)

MRC

For Help regarding configuration of the syngo MR system, please click **Help** at the top right of the Configuration platform refer to "Procedures MRC".

- After completed configuration click "Finish", re-boot the system as follows: Click "**Home**", confirm the re-boot message window with "**OK**", **please wait the system will need ~5 minutes before it will re-boot.**
- Reconnect the customer LAN connection

NOTE

Check / install Hotfixes and Bugfixes

Activate BIOS protection

- Insert BIOS CD "MR059 BIOS CD VA25E with PN: 100 49 286 P9901"
- Reboot the system
- Wait until: "HIPAA installed" is displayed
- Remove the BIOS CD

IFN (Image Filter New) Installation Work Steps

NOTE

The IFN is only available for newly delivered VA25A-based MR systems or Evolve Upgrades beginning with 2006 and up. In case of a software reinstallation, check if the following CD "syngo MR 2004A, Numaris/4 VA25A - IFN Rev. A, CD 1/1, Part no.: 10018855 K2219" is available on the system. If not, the installation of the IFN can be skipped for the software reinstallation

NOTE

For newly delivered VA25A-based MR systems or Evolve Upgrades beginning with 2006 and up, the IFN software replaces the image filter software from Context Vision. The filter function of the IFN (Image Filter New) is the same as that of the old one from Context Vision.

NOTE

The IFN can be installed only after the service pack 2 update.

NOTE

The IFN can be installed only before the service pack 3 update.

Hints

- On some systems, the following or a similar error message may be displayed occasionally:



Fig. 15:

Installing IFN on MRC and MRSC

ACTION	MRC	MRSC
- Check if SoftwareUpdate VA25A-SP2 is already installed. If not, please refer to the UI MR006/05/R or MR008/05/P. In the main syngo MR UI, you can verify if the SP2 is installed: - Help -> Info - Check if the "N4_VA25A_LATEST_20040724_P16" entry is visible under Siemens Software Packages or see logfile C:\MedCom\MriProduct\inst\SWHistory.txt	☐	☐
To install the IFN package, you have to be logged in as Administrator. If necessary, reboot and log in as Administrator (disable the autologin function in the local service UI).	☐	☐
Insert the CD "syngo MR 2004A, Numaris/4 VA25A - IFN Rev. A, CD 1/1, Part no.: 10018855 K2219" in the CD-ROM drive.	☐	☐
- Start Window Explorer: - Click START -> Programs -> Accessories -> Windows Explorer	☐	☐
Navigate to drive R:	☐	☐
Run program MR_Update.bat by double-clicking on it. NOTE: The program will first copy the software update package to the system. After copying is complete, installation of the package will start in a new CMD window (this will take approx. 20 - 30 min!) The original CMD window displays 'Press any key to continue ...' This can be done safely at any time. The CD will be ejected AUTOMATICALLY. Logfile entries such as: "Installation finished with ERRORLEVEL: 42" are indicating no errors!! The installation has finished when the following message is displayed in the CMD window entitled 'Software Update Installation - C:\sysmgmt\sd_store\MR_Update.bat INST': Finished software update at "actual date" "actual time", please reboot system!	☐	☐
Wait until please reboot system! is displayed.	☐	☐

ACTION	MRC	MRSC
Remove the CD "syngo MR 2004A, Numaris/4 VA25A - IFN Rev. A, CD 1/1, Part no.: 10018855 K2219" from the CD-ROM drive.	☐	☐
Open a cmd tool window.	☐	☐
Type in the following commands: SET TARGET=%Medhome%\MriProduct\ExamDbNls\en move /y %TARGET%\MrExamDb_EdbResen.dll %TAR- GET%\MrExamDb_EdbResen_old.dll copy /y %TARGET%\MrExamDb_EdbRes.dll %TAR- GET%\MrExamDb_EdbResen.dll	☐	☐
Reboot the complete system.	☐	☐

Verifying correct IFN installation

ACTION	MRC	MRSC
- In the main syngo MR UI, you can verify correct installation as follows: - Help -> Info - Check if the "N4_VA25A_LATEST_20040724_P25" entry is visible under Siemens Software Packages. If not, the IFN is not installed correctly.	☐	☐

CV License (only valid for systems delivered before 01/2006)

- If you have the option CV license filter, the low level key must be entered in the Service UI Utilities, Escape to OS with the command: **cvlicense** [license key]. For detailed information see ([Requested data for cvlicense generation / p. 141](#)).

What to do next?	
	Refer to the flowchart (Fig. 2 / p. 51) after finishing the configuration

MRSC

For Help regarding configuration of the syngo MR system, please click **Help** at the top right of the Configuration platform refer to "Procedures MRSC"

General

NOTE

An empty CD-R is used as medium for the procedure ([Export MRC configuration data to MRSC / p. 59](#)) and ([Import configuration data from CD-R to MRSC / p. 60](#))

(not a floppy as described in the On-line Help texts)

NOTE

SOFTWARE ERROR IMPORTANT!!!

When you perform Initial configuration

Database

Satellite Console

DON'T CLICK SAVE!

Refer to below workflow:

Workflow to avoid unwanted re-boot		
Local Host		
	System Type	
	Site Info	
	TCP/IP LAN	
	Users	
	MonitorType	
Database		
	Satellite Console	DON'T click SAVE just use > to navigate
Service		
	Licensing	Import your MRSC licenses
Navigate back to "Database"		

Workflow to avoid unwanted re-boot		
	Satellite Console	Perform configuration for your type of satellite console, save the system will re-boot automatically when you click save
MRSC shares DB with MRC: See (MRSC Sharing Database with MRC / p. 72)		or Local DB at MRSC: See (Completed configuration / p. 72)

MRSC Sharing Database with MRC

For Help regarding configuration of the syngo-MR system, please click Help at the top right of the Configuration platform refer to "Procedures MRSC"

Completed configuration

- After completed configuration click "Finish", re-boot the system as follows: Click **"Home"**, confirm the re-boot message window with **"OK"**, **please wait the system will need ~7 minutes before it will re-boot.**
- Reconnect the LAN connections

NOTE

Check / install Hotfixes and Bugfixes

- Make sure that software options with additional installation media (CD's or DVD's) are installed

Activate BIOS protection

- Insert BIOS CD "MR059 BIOS CD VA25E with PN: 100 49 286 P9901"
- Reboot the system
- Wait until: "HIPAA installed" is displayed, remove BIOS CD from the CD

IFN (Image Filter New) Installation Work Steps

NOTE

The IFN CD is only available for newly delivered VA25A-based MR systems beginning with 2006 and up. In case of a software reinstallation, check if the following CD

"syngo MR 2004A,

Numaris/4 VA25A - IFN Rev. A,

CD 1/1,

Part no.: 10018855 K2219"

NOTE

is available on the system. If not, the installation of the IFN can be skipped for the software reinstallation

NOTE

For newly delivered VA25A-based MR systems or Evolve upgrades beginning with 2006 and up, the IFN software replaces the image filter software from Context Vision. The filter function of the IFN (Image Filter New) is the same as that of the old one from Context Vision.

NOTE

The IFN can be installed only after the service pack 2 update.

NOTE

The IFN can be installed only before the service pack 3 update.

Hints

- On some systems, the following or a similar error message may be displayed occasionally:



Fig. 16:

Installing IFN on MRC and MRSC

ACTION	MRC	MRSC
- Check if SoftwareUpdate VA25A-SP2 is already installed. If not, please refer to the UI MR06-05-R or MR08-05-P. In the main syngo MR UI, you can verify if the SP2 is installed: - Help -> Info - Check if the "N4_VA25A_LATEST_20040724_P16" entry is visible under Siemens Software Packages or see logfile C:\MedCom\MriProduct\inst\SWHistory.txt	☐	☐
To install the IFN package, you have to be logged in as Administrator. If necessary, reboot and log in as Administrator (disable the autologin function in the local service UI).	☐	☐
Insert the CD "syngo MR 2004A, Numaris/4 VA25A - IFN Rev. A, CD 1/1, Part no.: 10018855 K2219" in the CD-ROM drive.	☐	☐
- Start Window Explorer: - Click START -> Programs -> Accessories -> Windows Explorer	☐	☐
Navigate to drive R:	☐	☐
Run program MR_Update.bat by double-clicking on it. NOTE: The program will first copy the software update package to the system. After copying is complete, installation of the package will start in a new CMD window (this will take approx. 20 - 30 min!) The original CMD window displays 'Press any key to continue ...' This can be done safely at any time. The CD will be ejected AUTOMATICALLY. Logfile entries such as: "Installation finished with ERRORLEVEL: 42" are indicating no errors!! The installation has finished when the following message is displayed in the CMD window entitled 'Software Update Installation - C:\sysmgmt\sd_store\MR_Update.bat INST': Finished software update at "actual date" "actual time", please reboot system!	☐	☐
Wait until please reboot system! is displayed.	☐	☐

ACTION	MRC	MRSC
Remove the CD "syngo MR 2004A, Numaris/4 VA25A - IFN Rev. A, CD 1/1, Part no.: 10018855 K2219" from the CD-ROM drive.	☐	☐
Open a cmd tool window.	☐	☐
Type in the following commands: SET TARGET=%Medhome%\MriProduct\ExamDbNls\en move /y %TARGET%\MrExamDb_EdbResen.dll %TAR- GET%\MrExamDb_EdbResen_old.dll copy /y %TARGET%\MrExamDb_EdbRes.dll %TAR- GET%\MrExamDb_EdbResen.dll	☐	☐
Reboot the complete system.	☐	☐

Verifying correct IFN installation

ACTION	MRC	MRSC
- In the main syngo MR UI, you can verify correct installation as follows: - Help -> Info - Check if the "N4_VA25A_LATEST_20040724_P25" entry is visible under Siemens Software Packages. If not, the IFN is not installed correctly.	☐	☐

CV License (only valid for systems delivered before 01/2006)

- If you have the option CV license filter, the low level key must be entered in the Service UI Utilities, Escape to OS with the command: **cvlicense** [license key]. For detailed information see ([Requested data for cvlicense generation / p. 141](#)).

What to do next?	
	Remember to create a backup if the configuration has changed (refer to (Save site specific data / p. 79))

Overview installing Imager SW VA25

General

This chapter describes how to install the Imager software.

NOTE

The MAGNETOM syngo-MR Imager is delivered with factory-installed software. The following Installation procedure is necessary only in case of replacement of the host.

NOTE

Before performing one of the following steps, refer to [\(MR-000.840.12 / Host / Imager\)](#) for troubleshooting information.

Installing Software at the Imager "R630"

Prerequisites

- Customer-CD 4/6 MR 2004A Rev. B
- The measurement network is connected

Check/Import the Imager license to the MRC

- Login as administrator (MRC)

Use Windows file Explorer to copy/paste the license file from the floppy disk to the following target on the MRC: **C:\MedCom\tftp**

- Rename the license file for the Imager to **license.dat**

Installation procedure

- Insert the CD "Customer-CD 4/6." in the CD-ROM drive
(Windows XP, syngo, and syngo MR is included at CD 4/6)
- Switch off / on the power to force a re-boot from the CD-ROM drive
- Please wait ~10 minutes, the CD-ROM drive opens the carriage automatically, remove the CD.
- The silent installation will need ~45 minutes additionally

Following table is only provided as additional information for sites with a monitor and keyboard connected at the imager:

Output	Note	
A pseudo-graphic surface appears with the warning that all data will be lost after pressing a key.		All data will be lost at the system disk of the Imager. (Installation time ~10 minutes)
The Imager will reboot		
The CD-ROM drive opens the carriage automatically for CD 4/6		Remove the CD from the CD-ROM drive
Automatic reboot Windows XP Professional Setup Please wait while Setup copies files to your hard disk	Do not interfere the automatic installation / configuration. Windows XP Professional Setup will follow Several re-boots will follow Estimated time for the automatic configuration ~ 45 minutes	

NOTE**The Imager license must be available at the MRC**

- Boot the MRC (Host).
- Start the MR application.
- Switch on power to the Imager

NOTE**The Imager will reboot twice following the installation.****The initialization of each raw-data disk takes approximately 5 minutes!****(No activity will be seen in the task manager)**

(To check the measurement network connectivity to the Imager enter the command **ping -t mrir** in a cmd tool at the MRC)

- The imager is ready when: on each image drive F:\ (at High-End extension additionally G:\, H:\, I:\) there is a file with the size of 4 GB placed with the name "rawexchange.dat".

The imager is ready to use, when this file exists on all the image drives with the above-mentioned size.

(Additional files at F: "MattauchBuffer1.dat" and G: "MattauchBuffer2.dat") (You can use the Remote Desktop Connection to access the Imager: Start, Programs, Accessories, Communications, Remote Desktop Connection. The computer name of the Imager is **mrir**), User: **MRSystem**.

What to do next?**Install all released Service Packs according Speedinfo Software Reinstallation Steps for VA25A****Make sure that the system is producing images.**

Save site specific data

Save the site-specific data each time you make changes to the site-specific configuration, tune-up, or sequences and protocols. In case of a hard disk crash, you will be able to restore a valid set of site-specific data.

Perform a backup of the following packages:

Numaris4

The Numaris4 package contains the complete system tune-up data, e.g. shim files, coil tune-up, measurement system-specific data

CustomerProtocols

The CustomerProtocols package contains the customer protocols.

SW-Settings

The SW-Settings02 package includes configuration data regarding e.g. remote nodes, DICOM nodes, Licensing, Customer UI...

- Table entries highlighted in gray are valid only in certain cases. Refer to the “Note”.

This part of the table is always valid.
This part of the table is only valid for special cases. Refer to the related “Note”.

NOTE	Backups to local disks can be performed only at the customer site, not via remote logon.
-------------	---

Prerequisites

You need a empty recordable compact disk medium.

NOTE	Only empty CD-R (recordable compact disk) media may be used to save site-specific data.
NOTE	If there is a HC Device "DSS" in the backup data set, the restore of the camera settings will fail.

Settings that are not part of the backup/restore

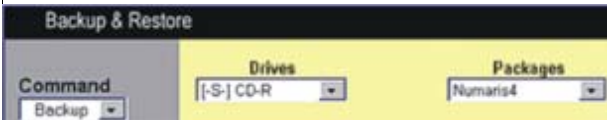
The following settings are not part of the site specific data. Make sure to note the settings in your configuration pages.

- Local Time zone
- **Local host, Users, "Enable Autologin"**
- Basic Worklist SCP Query Properties, Query Waiting Time
- **Printer installation, printer configuration**
- Static routes (See: [\(Check / add route to SIEMENS Remote Services router / p. 98\)](#))
- System language and weight
- **Settings for "Remote Service Access Control** (e.g., "Permanent Limited Access")

Backup (Save site specific data) on CD

- Login as **administrator** ↵.
- Start the service platform, select, **Backup & Restore**.

Output	Input	Note
 Fig. 17:	Select Backup in the Command pull-down menu	Select Backup
 Fig. 18:	Select the backup device in the Drives pull-down menu, e.g. CD-R	The CD-Burner has the drive letter S:
CustomerProtocols	Select Customer-Protocols in the Packages pull-down menu Click Go	Estimated backup time is ~ 1-12 minutes

Output	Input	Note
SW-Settings02	Select SW-Settings02 in the Packages pull-down menu Click Go	Estimated backup time is ~ 1 minutes
	Wait for message: "End Backup"	
 Fig. 19: Numaris 4	Select Numaris4 in the Packages pull-down menu Click Go	Estimated backup time is ~ 3-15 minutes
	Wait for message: "End Backup"	
SW-Settings02	Select MriCVLi-censeFile in the Packages pull-down menu Click Go	Estimated backup time is ~5 seconds
	Wait for message: "End Backup"	

Backup (Save site specific data) on USB device

NOTE

Please check the external media (e.g. USB-disk) first with the virus scanner of the service PC. Be sure that the external media used is free of viruses.

- Login as **administrator** ↴.
- Start the service platform, select, **Backup & Restore**.
- Start the "Internet Explorer"
- Connect the USB device.

Output	Input	Note
	Select the USB device in the Drives pull-down menu, e.g. F	E.g. drive letter F:
CustomerProtocols	Select Customer-Protocols in the Packages pull-down menu Click Go	Estimated backup time is ~ 1-8 minutes

Output	Input	Note
SW-Settings02	Select SW-Settings02 in the Packages pull-down menu Click Go	Estimated backup time is ~ 1 minutes
	Wait for message: "End Backup"	
SW-Settings02	Select Numaris4 in the Packages pull-down menu Click Go	Estimated backup time is ~ 2-10 minutes
	Wait for message: "End Backup"	
SW-Settings02	Select MriCVLicenseFile in the Packages pull-down menu Click Go	Estimated backup time is ~5 seconds
	Wait for message: "End Backup"	

Additional packages

Output	Input	Note
The following packages have to be backed up in certain cases only:		
Security-Settings Only valid for systems with HIPAA	Select Security-Settings in the Packages pull-down menu Click Go	Only to be performed at the MRC (and stand-alone MRSC's with a Local DB) . Security and user / group settings of HIPAA configuration

Output	Input	Note
Diagnostic Performed only in special service cases!	Select Diagnostic in the Packages pull-down menu Click Go	Files generated with MrSaveLog.bat Estimated back-up time is less than 1 minute
ExamDB Only valid for application specialists for special cases!	Select ExamDB in the Packages pull-down menu Click Go	Only to be performed at the MRC Estimated backup time 1-25 min, depending on the number of protocols

Restoring site specific data MR 2004A Rev. B on MR 2004A Rev. B

Restoring the site-specific data may be helpful in many situations. Please read the following information before beginning work.

- Troubleshooting: In case of corrupted site-specific data, we recommend saving the actual state of the site-specific data to another medium before restoring an old save set. This will save time if the problem is resolved by reinstalling the software. The actual tune-up data may save a great deal of work, especially if the customer recently optimized his protocols.
- Disk crash: After installing the software, a set of site-specific data may prove to be helpful.

NOTE

Backups to local disks can be performed only at the customer site, not via remote logon.

Prerequisites.

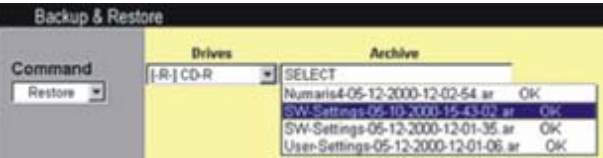
NOTE

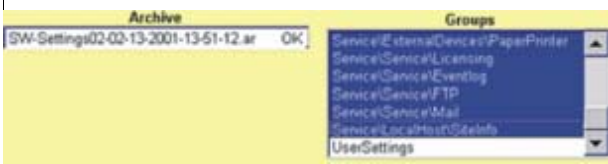
Ensure that the save set contains data valid for the current software version

Restore site-specific from CD.

(In case of initial configuration, a security warning is displayed. Confirm with "Yes" and use **Alt-Tab** to toggle to the Local Service Platform)

- Login as user: **administrator**
- Insert Backup CD-R, Select, **Backup & Restore** in the service platform

Output	Input	Note
 Fig. 20:	Select Restore in the pull-down menu Command Select device, e.g. CD-R the pull-down menu Drives	The CD-ROM will get the drive letter R:
 Fig. 21:	Select SW-Settings02 in the Archive pull-down menu	

Output	Input	Note
 Fig. 22:	Select all groups except " User-Settings and without SecConfig)" in the Groups pull-down menu	Due to a software bug, the User-Settings group must be restored separately after a system reboot
General use of mouse and key combinations To select all groups at once, click the group at the top, press and hold the Shift key, and click the group at the bottom. All groups are highlighted in blue. (If you press and hold the Ctrl key, any combination of groups can be selected by clicking the mouse)		
	Click Go	

- Wait until the message "***** End Restore:" is displayed

Output	Input	Note
	Select MriCVLicenseFile in the Archive pull-down menu	Low level license of Context Vision filter
	Click Go	

- Wait until the message "***** End Restore:" is displayed

Output	Input	Note
	Select Numaris4 in the Archive pull-down menu	Ignore error message "License file ...Error"
	Select all Groups	Confirm the popup window: 'Errors in RestoreHook' with OK
	Click Go	

- Wait until the message "***** End Restore:" is displayed
- Wait until the message "***** End Restore:" is displayed

Output	Input	Note
	Select CustomerProtocols in the Archive pull-down menu	Only restore the package "CustomerProtocols" if it's more recent than the package "Numaris4" (The package Numaris4 contains the customer protocols as-well.)
	Select all Groups	
	Click Go	

- Wait until the message "***** End Restore:" is displayed
- Click "**Home**" wait for the report generation to finish. Confirm pop-up with "**Cancel**"

Restore site-specific from external USB device

NOTE

Please check the external media (e.g. USB-disk) first with the virus scanner of the service PC. Be sure that the external media used is free of viruses.

Restoring the site-specific data may be helpful in many situations. Please read the following information before beginning with the work.

In case of corrupted site-specific data, we recommend saving the actual status of the site-specific data to another medium before restoring an old save set. This will save time if the problem is resolved by reinstalling the software. The actual tune-up data may save a great deal of work, especially if the customer optimized his protocols recently.

- Login as **administrator** ↵.
- Start the "Internet Explorer"
- Login to the service platform, select **Backup & Restore**.
- Connect the USB device.
- Select **Restore** in the "Command" drop-down menu.
- Select the source USB-drive (drive letter e.g. **F:**) in the "Drives" menu.
- Select your packages to be restored xxx_xxx_xxxx_xxxx_ok.ar in the "Archive" menu.
- In the menu "Groups", select the necessary packages to be restored.
- Click **Go** to start the "Backup", check for errors in the log window after restore has been completed.

NOTE

To ensure the Administrator account password is correct following restoration of the dynamic data, it is recommended to set the default Administrator password in the **Service GUI > Configuration > Local Host > Users**

Settings that must be manually activated after restore

Following configuration data needs manual actions to be activated after a restore has been performed.

- Local Time zone
- In case of host replacement the local TCP IP settings are not restored.
- **Local host, Users, "Enable Autologin"**
- Basic Worklist SCP Query Properties, Query Waiting Time
- **Printer installation, printer configuration**
- Static routes (See: [\(Check / add route to SIEMENS Remote Services router / p. 98\)](#))
- System language and weight
- **Settings for "Remote Service Access Control** (e.g., "Permanent Limited Access")

NOTE

In case the error message "C:\MedCom\MriSitedata\ContextVision\siemensIRU-PCNT-init-> access denied" is displayed during restore of N4 data, re-enter the cvlicense as described in (Context Vision / p. 141).

- Activate the measurement system configuration data as follows: Start the local Service platform, select **Configuration**, and click **Save** in each platform listed under **measurement** (the changes will be effective after the next reboot).
- Reboot the system from the **start** menu, select **Start -> Shut Down -> select Restart** and press **OK**.
(in case of SW reinstallation, the reboot takes ~ 20 additional min.)
- Login as **meduser**
- Select **"SW-Settings02"** in the Archive pull-down menu. Select **"User-Settings"** in the Groups menu.
- Click **"Go"** to start the restore. Read the messages in the log window.
- Wait until the message "***** End Restore:" is displayed.
- Click **"Home"** wait for the report generation to finish.
- Confirm pop-up window with **"Cancel"** and then reboot system manually after closing the Service Platform - Reboot the system from the **start** menu, select **Start -> Shut Down -> select Shutdown** and press **OK**.
- **Re-boot the system** (Connect the customer LAN if disconnected)
- After the reboot, perform a spot check regarding image quality. Also, verify that all customer options are listed in the system imprint.

NOTE

The HIPAA configuration is only to be restored at the MRC (and standalone MRSC's with a Local DB).

Output	Input	Note
Before restore of security settings a reboot of the Host is necessary.		
		Login as meduser
Security-Settings Only valid for systems with HIPAA	Select Security-Settings in the Packages pull-down menu Click Go	Only to be performed at the MRC (and stand-alone MRSC's with a Local DB) . MRSC that shares the DB with MRC also shares the security settings with the MRC

General

For systems that include a satellite console, **configure the main console first.**

NOTE

You need your notes from the configuration pages.



The configuration pages contain important information regarding the default configuration/settings for syngo MR systems refer to ([Configuration pages for MRC / p. 20](#)), for MRSC refer to ([Configuration pages MRSC / p. 40](#)).

NOTE

MRC/MRSC are used as abbreviations for the main/satellite console, mrc/mrsc are used as host names.

NOTE

Enabling of security settings (HIPAA) is to be performed by quali-fied personnel with HIPAA training.

NOTE

The use of non-released combinations of hardcopy cameras and BU products is prohibited and does not meet the requirements defined. For detailed information regarding the CAMERA configuration, refer to technical documentation "Camera Info Numaris 4 VA25A"

NOTE

The "SIEMENS Remote Services TM" will only support router connections for NEWLY delivered systems. See <http://cs.med.siemens.de/> "For Service, Siemens Remote Services" for the most recent information regarding router configuration.

NOTE

If you have the option "Image Filter Software"(Context vision) an additional license must be enabled.

Refer to ([Context Vision / p. 141](#)).

NOTE

For known bug, please refer to "Known bugs and hints for Service"

NOTE

It's recommended not to configure WINS or DNS for syngo MR systems. (Netstat -a returns putpunt ~1h when WINS and DNS are configured)

Prerequisites

- The configuration pages contain the information/notes required for your system.
- The VA25A software is installed.

Access Service Software

- Simultaneously press the **Ctrl - Alt - Del** keys to open the login window.
- Login as user **meduser**.
- Click **Options** and select **LocalService**.

Output	Input
Security warning	Click Yes
 The image shows the Siemens Service Software Authentication window. At the top is the Siemens logo. Below it, the text 'Authentication' is in a black box. The main title is 'Service Software'. Below that, it says 'Erlangen', 'Harmony MR 11002', and 'NUMARIS/4'. A prompt 'Please enter password' is followed by a password field containing '9CA8DAF5B76323' and an 'OK' button. At the bottom, there is a checkbox labeled 'Set as Default'.	Remember to click Set as Default Enter the variable key Click OK
Fig. 23:	Click Next

NOTE

For detailed information regarding the configuration of syngo MR systems, refer to "Help" in the Service Software Platform.

Changes in the configuration

- If you have changed the configuration, re-boot the system as follows: Click "Home", confirm the re-boot message window with "OK"

What to do next?	
	Remember to create a backup if the configuration has changed (refer to (Save site specific data / p. 79)).

Parallel Printer Installation

NOTE

System with two consoles (MRC and MRSC), local printer connected via the parallel port is only supported at the MRSC.

- Login as user Administrator

Installing local parallel printer

- Select Start, **Settings, Control Panel**.
- Click **Printers and Other Hardware**
- click **Printers and Faxes**
- Click the **Add Printer** icon.
- Click **Next** in the Add printer Wizard "Welcome to Add Printer Wizard"
- Select radio button **Local printer attached to this computer**, click **Next**
- Ignore that the Add Printer Wizard "Searching..." terminates with the message "The wizard was unable to detect any Plug and Play printers..."
- **Select LPT1 ... in the pull-down menu "Use the following port:", click Next**
- Select your printer in the Add printer Wizard "Install Printer Software", click **Next**
- Select radio button **Keep existing driver (recommended)**, click **Next**
- Name the printer in the Add printer Wizard "Name Your Printer", click **Next**
- Give a logical Name for the printer in the Add printer Wizard "Printer Sharing", click **Next**
- Select radio button **Yes** in the Add printer Wizard "Print Test Page", click **Next**
- Click the **Finish** in the printer Wizard "Completing the Add Printer Wizard".
- Click **OK** in the pop-up "A test page is now being sent to the printer ..."

Enabling MRC to access printer share of MRSC

NOTE

If the actual console is MRSC, you must grant access to the MRC in the IP security filter.

- Keep login as user Administrator (at the MRSC)
- Open a cmd window
- enter the command **ipseccmd -f (0+146.254.111.222) exchange 146.254.111.222** in this example with the IP address of your MRC

Network Printer

The following procedure must be performed at each MR console that will print to the shared network printer.

- Login as user Administrator

Installing network printer

- Select Start, **Settings, Control Panel**.
- Click **Printers and Other Hardware**
- click **Printers and Faxes**
- Click the **Add Printer** icon.
- Click **Next** in the Add printer Wizard "Welcome to Add Printer Wizard"
- Select radio button **A network printer, or a printer attached to another computer**, click **Next**
- **Select radio button Connect to this printer (or to browse for a printer ... , enter the IP address of the MRSC e.g. \\146.254.111.221 in field "Name:" (In the Add printer Wizard "Specify a Printer"), click Next**
- Click **Yes** in the pop-up "Connect to Printer" You are about to connect to a printer on xx.xx.xx.xx ...
- Click the **Finish** in the printer Wizard "Completing the Add Printer Wizard".

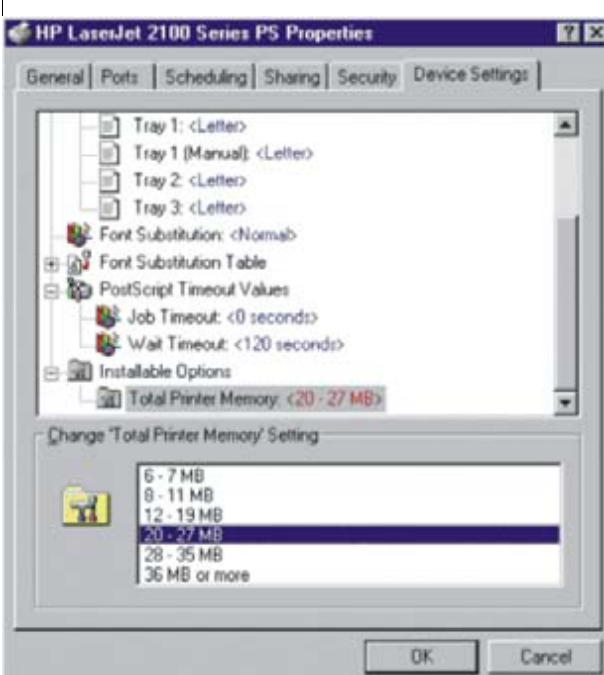
General Printer Properties

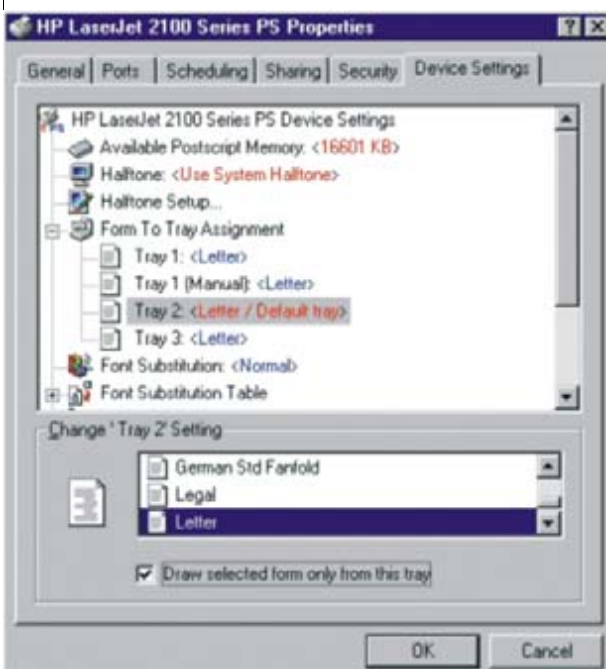
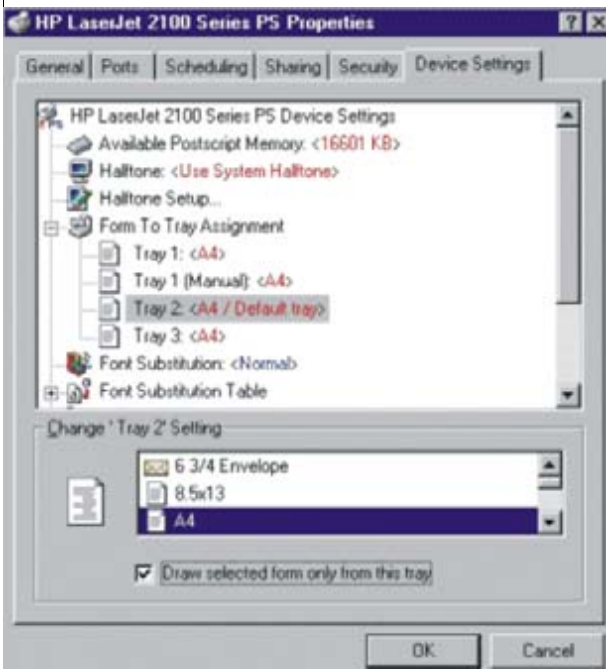
Each MR console that will be used to print has to be configured as follows to optimize print quality.

NOTE	The printer should have at least 24MB RAM.
-------------	---

Output	Input	Note
	Click Start and select Settings, Printers	Click START in the lower left corner
	Right-click the HP printer icon and select Properties	

Fig. 24:

Output	Input	Note
 Fig. 25:	Click the Device Settings tab Navigate to Total Printer Memory Adjust the "Total Printer Memory Setting"	A minimum of 24 MB are recommended
 Fig. 26:	Check that half tone is set, select "Use System Halftone" radio button	

Output	Input	Note
 Fig. 27: Paper size: US Letter	Navigate to "Form To Tray Assignment" Click Tray2: (Select used paper size for the customer) Click the Draw selected form only from this tray check box	Only valid for printers using the "Letter" paper size
 Fig. 28: Paper size: A4	Navigate to "Form To Tray Assignment" Set each tray to A4 Click Tray2: Click the check-box Draw selected form only from this tray check box	Only valid for printers using the "A4" paper size
	Click OK	

Output	Input	Note
	Right-click the HP printer icon and select Document Defaults	
	Click the Advanced tab Letter or A4 For A4 paper size, select A4 For "Letter", no changes are necessary Navigate to "Graphic" Select "Resolution" Set the resolution to 600 dpi	
		Click OK
Reboot the system		

Fig. 29: Document Defaults

Fig. 30: Paper size & Resolution

- Remain logged in or log in as user "administrator"
- Log in as user "meduser".
- Start the local service platform, select **paper printer** under "external devices" and set the paper printer media formats available for your printer.
- Print a test image from the application to verify the print quality.

General

Only IE 6 has been tested and released for remote service as service browser at syngo MR systems running with NUMARIS MR 2004A Rev. B.

The project manager is responsible to plan for a router connection to the Siemens Remote Diagnostic server. Regarding questions to your pre configured router from Siemens please contact (mailto:srs_final@med.siemens.de, mailto:srs_final@med.siemens.de) or your remote diagnostic responsible.

Prerequisites

- Router is pre-configured and known as route to SRS by the MRC (MRSC if available)
- System is configured at the SRS server

System Management

System management is a new remote service for syngo MR systems, it provides following functionality: proactive event reporting and statistical data transfer. In the future it will be used for SW distribution e.g hotfix / buxfix and anti virus patterns.

Pro Active Service

With MR 2004A Rev. B we provide Autoreports for proactive event triggering, that can be distributed on demand.

Quality Assurance DB (QDB)

The QA data will be transferred to the server at regular intervals for statistic evaluation and trend analysis, the data is collected into the MR-quality data-base.

Helium level and magnet shield temperature data is transferred to the server at regular intervals. This function offers a quick overview of the current He levels and magnet shield temperatures for the reporting sites.

Under following URL you can see data that has been sent by your system: <http://td-quality.med.siemens.de/>

Event DB (EDB)

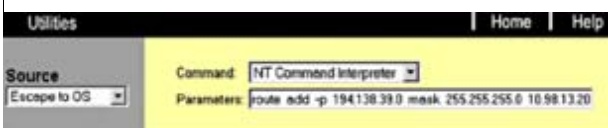
The EDB provides reports regarding events that are stored in the system event log.

System Utilization

Utilization Management is an new service offering and for leagal aspects the customer has to sign a contract in advance. Please involve your service organization before you activate System Utilization.

Check / add route to SIEMENS Remote Services router

As the most optimal solution, the customer network administrator has configured the route from the default gateway to the SIEMENS SRS server for your continent. In the case of a VPN connection being used for SRS, the customer network administrator must configure the route from the default gateway to the SIEMENS SRS server.

	Output	Input
 The screenshot shows the 'Route print' command output in the NT Command Interpreter. The output is divided into two sections: 'Interface List' and 'Active Routes:'. The 'Interface List' shows details for several network interfaces. The 'Active Routes:' section displays a table of active routes with columns for Network Destination, Netmask, Gateway, Interface, and Metric.	Start Local Service Select Utilities Select Escape to OS in the pull-down menu "Source" Enter the command in the Parameters field route print ↵ Check if the router for rdiag is listed In case the route is not listed to your SRS subnet you must add a static route.	
 The screenshot shows the 'Route add' command being entered into the Parameters field of the NT Command Interpreter. The command is 'route add -p 194.138.39.0 mask 255.255.255.0 194.138.39.0'.	The below example means use the customer router / gateway with IP: xx.xx.xx.xx for the SRS network as destination. E.g. (LUX09505 194.138.39.18) 194.138.39.0 route to the whole SRS subnet, and -p makes the entry persistent)	
(Replace the xx.xx.xx.xx with your customer router / gateway IP address !!!)		
Europe, Africa:	route add -p 194.138.39.0 mask 255.255.255.0 xx.xx.xx.xx ↵	
North-/ South America:	route add -p 129.73.116.0 mask 255.255.255.0 xx.xx.xx.xx ↵	
Asia:	route add -p 194.138.243.0 mask 255.255.255.0 xx.xx.xx.xx ↵	

Check to ensure that configured SIEMENS Remote Server is added to the IP security list, click Reports; select IP Security below "SiteInfo". If the System Management server is missing, refer to [\(Adding host to IP security list / p. 139\)](#).

Only in case you have a wrong route, use the following command to delete: route delete 1x.x.x.x ↵	(To delete a wrong route from the routing table: enter IP address of the routing network)
--	--

Configure MAIL and FTP target (router systems)

Output	Input	Note
 Fig. 33: Mail (name and IP of your SRS)	Enter hostname and IP address of your Access Server see below table: Select radio-button LAN Click Save Click OK in the pop-up window	

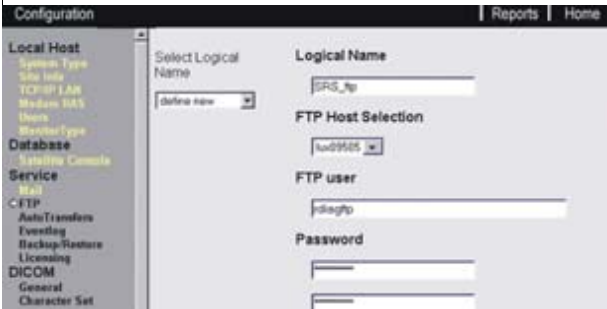
Country	Rdiag server name	IP of Access Server
Europe, Africa	LUX09505	194.138.39.18
North- / South America	SRSACC1	129.73.116.92
Asia	sgpt806x	194.138.243.178

For MR 2004A Rev. B it will be possible to import proactive templates from the Siemens Remote Server for syngo MR systems. You must use the logical names **SRS_mail** and **SRS_ftp** for the mail and ftp targets in order to be able use the preconfigured templates

Output	Input / Note
 Fig. 34: Logical name for mail target SRS_mail	Define mail target Logical Name SRS_mail Select your Access Server in the pull-down menu "Mail Server Host Selection" e-mail address, refer to the table below Click Save Click OK in the pop-up window

If you don't have a e-mail address for your USC, use one of the following address:

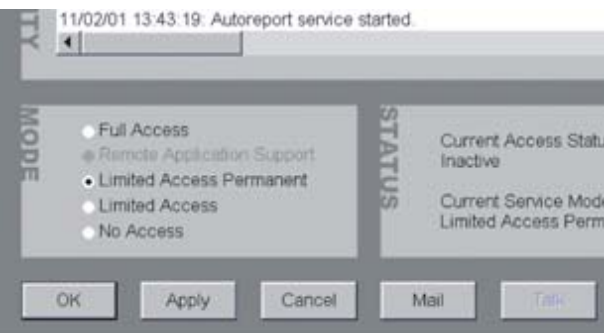
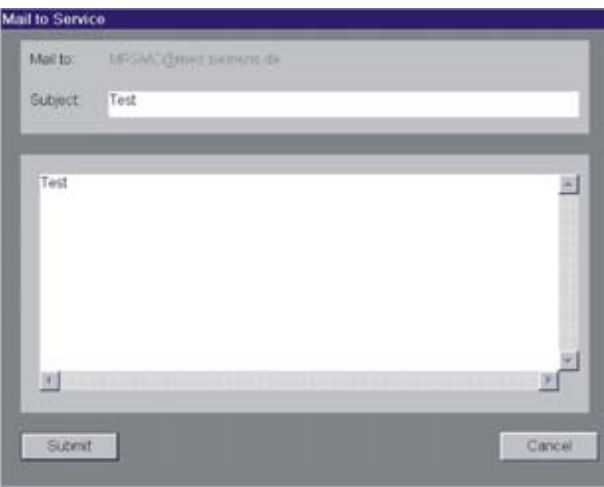
Access Server	e-mail address
Europe, Africa, and Asia	MRSMC.med@siemens.com
North- / South America	MRUPTIME@sms.siemens.com

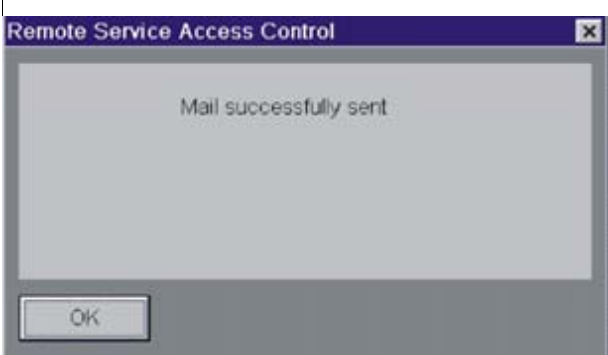
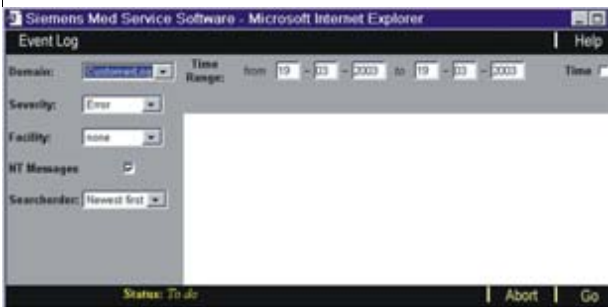
Output	Input	Note
	Skip the first window with >	The host target of mail is used for ftp as well.
 Fig. 35: Logical name for ftp target SRS_ftp	Define ftp target Logical Name SRS_ftp Select your Access Server in the pull-down menu "Ftp Host Selection" FTP user rdiagftp enter the password confirm the password Click Save Click OK in the pop-up window	
 Fig. 36: Activate automatic Transfers	If the customer agrees on remote monitoring Set the check-box: activate automatic Transfers	
For questions regarding the password for the rdiagftp account, please contact your remote diagnostic responsible / USC.		

Testing connectivity from system to Rdiag server

Send a test mail

- Check at the bottom right that no remote service is active

Output	Input / Note
 Fig. 37: Remote Service Access Control	Select Options, Remote Service Click Mail
 Fig. 38: Mail to Service	Enter subject e.g. SRS test MRxxxxx Text: Please confirm reception of this mail to your.name@siemens.com or under phone number xxx Click Submit

Output	Input / Note
 <i>Fig. 39: Mail successfully sent</i>	Click OK
 <i>Fig. 40: Mail successfully sent</i>	Select Options , Event Log Select Domain: Customer Log

If following error message is displayed

E 19-Mar-03 11:18:35 16 TransferMgr n/a MAIL error.

Receiver: "MRxxx@xxx.siemens.xx"

Filename: "C:\MedCom\service\extract\MRxxxxxx20030319111756Mail.txt".

the customer router is not properly configured for outgoing traffic, please contact the person responsible for the router configuration and the local network administrator.

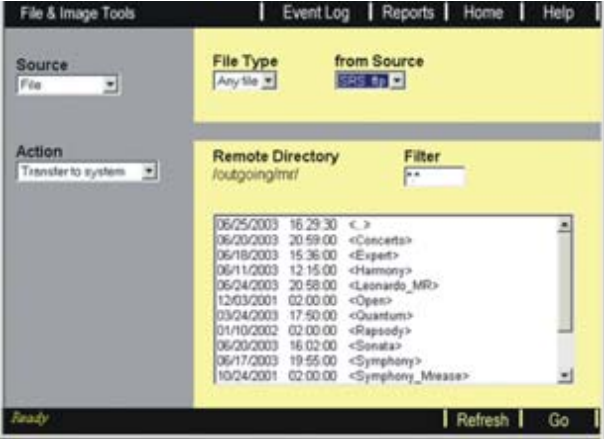
Check to ensure that the configured SIEMENS Remote Server is added to the IP security list, click **Reports**; select **IP Security** below "**SiteInfo**". If the System Management server is missing, refer to [\(Adding host to IP security list / p. 139\)](#).

Test FTP login to the Remote Server

Prerequisite

The test [\(Send a test mail / p. 102\)](#) was successful

- Start the Local Service platform, select File & Image Tools

Output	Input / Note
 Fig. 41: File & Image Tools	Source: select File Action: select Transfer to system File Type: select Any file from Source: select SRS_ftp
 Fig. 42: Listing "Remote Directory"	If a listing of folders is displayed in the frame "Remote Directory" the access to the FTP account of the Remote Server was successful. Click Home Close the Service Platform

- If the test fails, re-enter the password for account "**rdiagftp**" perform the test once more if the test fails again, please contact your Remote Server responsible / USC.

AutoReport configuration

It is important to transfer data from the site to your USC/TSC/HSC on a regular basis. Quality Assurance and Tune-Up measurements logs, which will be collected in a central database. This allows various kinds of statistic evaluation. The AutoReport function takes over the tedious job of regularly checking for new data and transmitting them manually.

NOTE

Be aware, that AutoReport can only work correctly, if the following terms are fulfilled:

LAN available ROUTER / Modem is connected and powered on

Your site has been entered into the RDIAG server

Target definition is correct. (Dial mode, escape number, phone

number and password (obtained from your USC/TSC/HSC)

AutoReport is defined and enabled.

The router is required free-of-charge to our customers during the warranty period. It also extends beyond that period if a service contract is concluded.

The planning of the "Remote Diagnostics" environment is performed by the project manager. For detailed information, please refer to the ("[Planning Guide](http://cs.med.siemens.de/cms/productinfo/default.asp?startfolder=801281ea3099989879a6fa4ba28d2b0c)", <http://cs.med.siemens.de/cms/productinfo/default.asp?startfolder=801281ea3099989879a6fa4ba28d2b0c>).

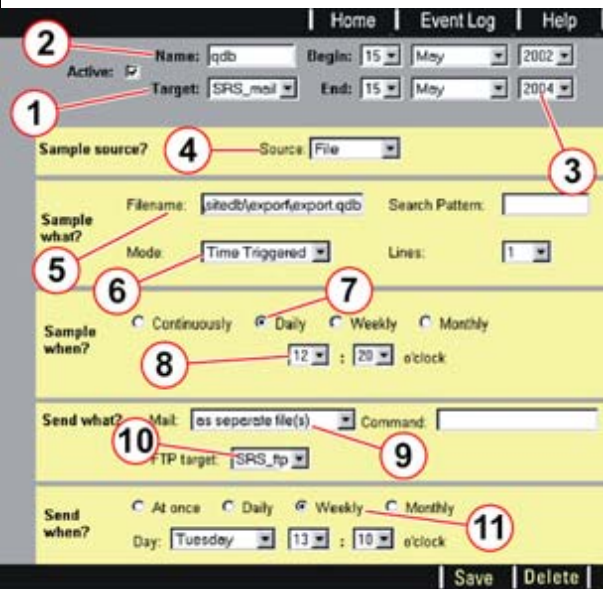
Configuring QDB report

- **Location of the QDB file: C:\MedCom\service\html\seso\sitedb\export\export.xml**
- **Sample When: Choose a time, when the system is usually in operation!**

The system saves the "AutoReport" configuration, then it sends a test-mail within 5 minutes after completed AR configuration to the mail target destination.

Set-up the Autoreport frequency to daily, to report magnet monitoring data.

- Start the service platform

Output	Input	Note
 Fig. 43:	Click at Auto Report Select TEMPLATE in pull-down menu REPORT	
 Fig. 44: Configure qdb auto report Pos. 1 Target Pos. 2 Name Pos. 3 End time Pos. 4 Source Pos. 5 Filename Pos. 6 Mode Pos. 7 Sample when Pos. 8 Sample when: time Pos. 9 Mail Pos. 10 FTP target Pos. 11 Send when	1. Target: Select "SRS_mail". 2. Name: e.g. qdb 3. End time e.g. 2011 4. Sample source: Select File 5. Filename: Enter location of the "qdb" file, see above path in "Location of the QDB file". 6. Sample what: "Mode:" Select Time Triggered 7. Sample when: Click radio button Daily 8. Sample when: time, see above "Sample When". 9. Send what?, "Mail:" Select as separate file(s) without mails 10. Send what?, "FTP target:" "SRS ftp" 11. Click radio button At once	
When all of it is configured correctly, click the Save button.		
The system will send a test-mail automatically after saving the AutoReport For results refer to (Send a test mail / p. 102).		

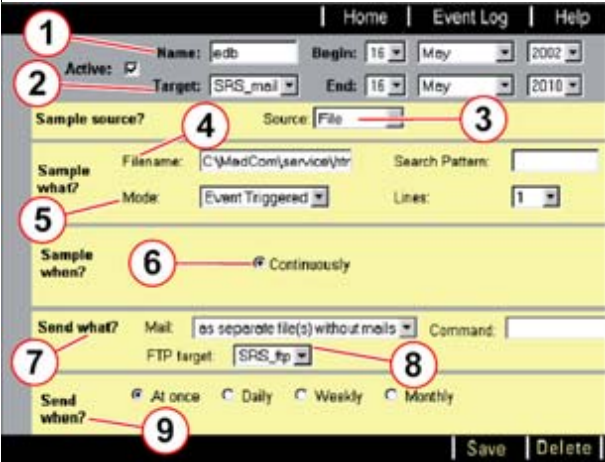
Configuring EDB report

- **Location of the EDB file: C:\MedCom\service\html\seso\sitedb\export*.edb**

The system saves the "AutoReport" configuration, then it sends a test-mail within 5 minutes after completed AR configuration to the mail target destination.

Set-up the Autoreport frequency to once a week.

- Start the service platform

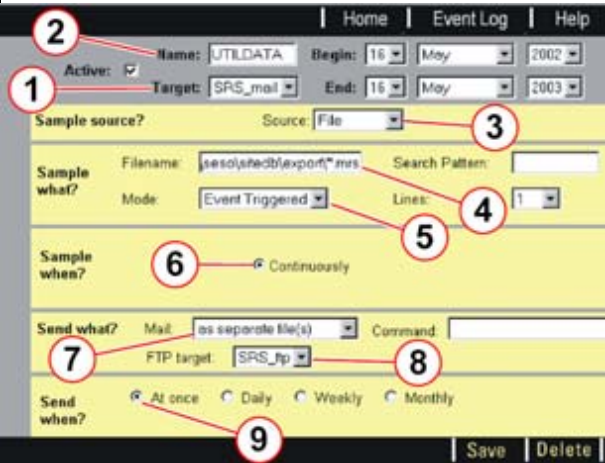
Output	Input	Note
REPORT TEMPLATE Fig. 45:	Click at Auto Report Select TEMPLATE in pull-down menu REPORT	
 Fig. 46: Configure edb auto report Pos. 1 Name Pos. 2 Target Pos. 3 Source Pos. 4 Filename Pos. 5 Mode Pos. 6 Sample when Pos. 7 Send what Pos. 8 FTP target Pos. 9 Send when	1. Target: Select "SRS_mail". 2. Name: e.g. edb. 3. Sample source: Select File 4. Filename: Enter location of the "mrs" files, see above path in "Location of the mrs file". 5. Sample what: "Mode:" Select Event Triggered 6. Sample when: Click radio button Continuously 7. Send what?, "Mail:" Select as separate file(s) without mails 8. Send what?, "FTP target:" Select SRSftp 9. Click radio button At once	
When all of it is configured correctly, click the Save button.		
The system will send a test-mail automatically after saving the AutoReport. For results refer to (Send a test mail / p. 102).		

Configure System Utilization reporting

- **Target:** Choose your responsible USC/TSC/HSC for support regarding rdiag server.
- **Location of sysutil files:** C:\MedCom\service\html\seso\sitedb\export*.zip

The system saves the "AutoReport" configuration, then it sends a test-mail within 5 minutes after completed AR configuration to the mail target destination.

- Start the service platform

Output	Input	Note
 Fig. 47:	Click at Auto Report Select TEMPLATE in pull-down menu REPORT	
 Fig. 48: <i>Configure System Utilization reporting</i> Pos. 1 Target Pos. 2 Name Pos. 3 Source Pos. 4 Filename Pos. 5 Mode Pos. 6 Sample when Pos. 7 Mail Pos. 8 FTP target Pos. 9 Send when	1. Target: Select rdiag host, see above "Target". 2. Name: UTILDATA 3. Sample source: Select File 4. Filename: Enter location of the "sysutil files" files, see above path in "Location of the file". 5. Sample what: "Mode:" Select Event Triggered 6. Sample when: Click radio button Continuously 7. Send what?, "Mail:" Select as separate file(s) without mails 8. Send what?, "FTP target:" Select SRS_ftp 9. Click radio button At once	
When all is correctly configured, click the Save button.		
The system will send a test-mail automatically after saving the AutoReport.		
For results refer to (Send a test mail / p. 102)		

System Utilization

NOTE

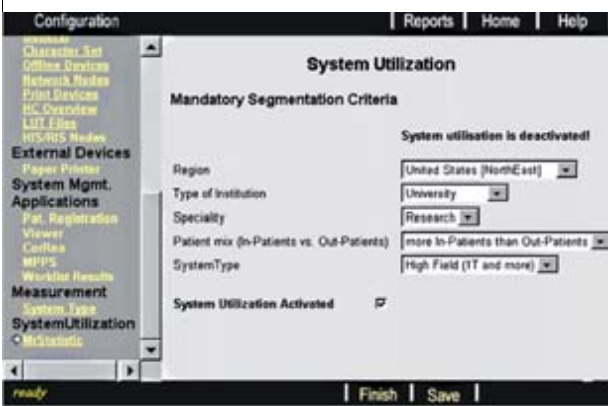
To meet legal aspects, the customer has to sign a contract in advance. Please involve your service organization before you activate System Utilization.

Autoreport for SYSUTIL must be configured before activating System Utilization.

NOTE

Only activate System Utilization for systems where the customer agrees on "Utilization Management"

The system utilization provides log and report data regarding the usage of syngo MR applications

Output	Input / Note
	Prerequisite to activate System Utilization is that activate automatic Transfers is checked in Configuration, AutoTransfers.
	Start Configuration in the Service UI, navigate to System Utilization
 <i>Fig. 49: Sytem Utilization Activated</i>	Use the pull-down menu to select region and site profile. If the customer has signed the contract on Utilization Management click the checkbox: System Utilization Activated Click Save to save changes

System Management Installation

System management is a new remote service for syngo MR systems, it provides following functionality: proactive event reporting and statistical data transfer. In the future it will be used for SW distribution e.g hotfix / buxfix and anti virus patterns.

The installation should be performed at the **MRC** and **MRSC** if available.

Prerequisites

Download the CA test tools from the **CS Intranet**: <http://cs.med.siemens.de>, **For Service, MR, Tools** download the **CA_Com_Check_Tool.exe** to your local service PC.

- Insert a blank floppy into the floppy drive
- Double-click the downloaded file to extract the file by clicking "Unzip". Check if 9 files are successfully extracted (popup message). Click Close to continue.
- Remove the floppy and label it "CA Com Test"

Enable IP security for the CA Unicenter servers in your region

Region: Europe and Africa	server name	IP address
Event Management Server	FTHW9MTA	194.138.39.20
Asset Management Server	FTHW9MUA	194.138.39.21
Software Distribution Server	FTHW9MVA	194.138.39.22

Region: North- / South America and Asia	server name	IP address
Event Management Server	SRSSQL01	129.73.116.90
Asset Management Server	SRSSQL02	129.73.116.93
Software Distribution Server	SRSSQL03	129.73.116.94

Enable the system management servers for your region to the IP security list, use following procedure:

- At the "MRC", simultaneously press **Ctrl - Esc**. Select **Start, Advanced User**; enter the advanced user password; click **OK**
- Simultaneously press **Ctrl - Esc**, click **Start, Run cmd** enter the command: **runas /user:administrator cmd** (enter the administrator password) In the newly started cmd tool you have administrator rights.

for Region: Europe and Africa type following commands in the cmd tool:
ipseccmd -f (0+194.138.39.20) ↵
ipseccmd -f (0+194.138.39.21) ↵
ipseccmd -f (0+194.138.39.22) ↵

for Region: North- / South America and Asia type following commands in the cmd tool:

```
ipseccmd -f (0+129.73.116.90) ↵
ipseccmd -f (0+129.73.116.93) ↵
ipseccmd -f (0+129.73.116.94) ↵
```

Routing

As the optimal solution, the customer network administrator has configured the route from the default gateway to the SIEMENS server for your continent.

- Skip this step if the SRS router is configured as default gateway. Do NOT add permanent routes if there is no entry under "Persistent Routes" of the logfile `routing.log` under `C:\Medcom\MrisiteData`)

Logfile example **without** any **permanent routes**:

```
.....
Persistent Routes:
    None
```

Logfile example **with permanent routes**:

```
.....
Persistent Routes:
    Network Address      Netmask                Gateway Address        Metric
      194.138.39.0       255.255.255.255       xxx.xxx.xxx.xxx        1
or
      129.73.116.0       255.255.255.255       xxx.xxx.xxx.xxx        1
or
      194.138.243.0      255.255.255.255       xxx.xxx.xxx.xxx        1
```


The below **example** means **use the customer router / gateway with IP : xx.xx.xx.xx** for the SRS network as destination. E.g. (LUX09505 **194.138.39.18**)
194.138.39.0 route to the whole SRS subnet, and **-p** makes the entry persistent)

Add permanent routes for system in region **Europe and Africa** to the **CA Server**:

```
route add -p 194.138.39.0 mask 255.255.255.0 xx.xx.xx.xx ↵
```

Add permanent routes for system in region **North- / South America and Asia** to the **CA Server**:

```
route add -p 129.73.116.0 mask 255.255.255.0 xx.xx.xx.xx ↵
```

(The Asia server "sgpt806x" automatically routes the IP addresses: 129.7.116.90, 129.7.116.93, and 129.7.116.94) via a VPN tunnel to the CA server)

Check connectivity to the System Management server

- Simultaneously press **Ctrl - Esc**, click **Start**, **Run cmd**

or use the following script from the floppy "**CA Com Test**" by entering:

for **Europe, Africa**

```
a: \CheckCA_Europe.BAT
```

for **North America, South America, Asia**

```
a: \CheckCA_US.BAT
```

NOTE

Important: Only continue the installation /configuration if the connectivity to the System Management server is successful.

Required ports for System Management (CA Unicenter) Communication

Provide following information to the customer network administrator **if the connectivity test failed** to the System Management server or contact your SRS responsible.

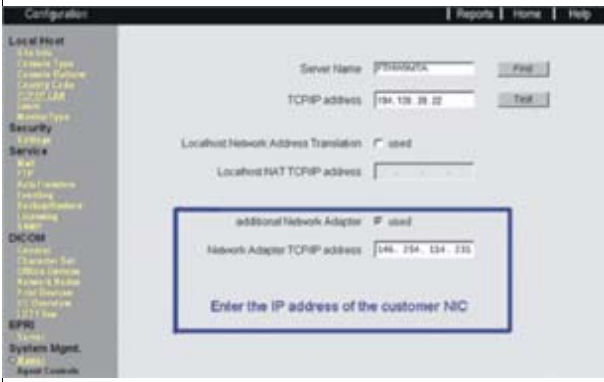
Protocol	Port	Description
TCP	4105	CA Messaging (TCP)
UDP	4104	CA Messaging (UDP)
TCP	8222-8235	Data Transport Agent (DTA)
TCP	7001	Common Communication Interface (CCI)
TCP	1721	CCI Remote Daemon

Protocol	Port	Description
TCP	4718-4724	Software Delivery (TCP)
UDP	4725	SD Job Check (UDP)

Master

Region: Europe and Africa	server name	IP address
Software Distribution Server	FTHW9MVA	194.138.39.22

Region: North- / South America and Asia	server name	IP address
Software Distribution Server	SRSSQL03	129.73.116.94

Output	Input / Note
 Fig. 50: CA Master	In this mask, you configure the Server Name / IP addresses of the System Management Master for your Region-see above table. Click the checkbox "additional Network Adapter if console type is MRC", enter the TCP/IP address of the network adapter used for SRS (Siemens Remote Services) in the field "Network Adapter TCP/IP address". This is usually the hospital network address of the host.
If NAT (Network Address Translation) is used, click checkbox "Local host Network Address Translation". As "Local NAT TCP/IP address", enter the NAT TCP/IP address assigned in the SRS server system configuration. Ask the local USC for this information.	

Output	Input / Note
Fig. 51: CA Functions	In this mask, you may stop and start the System Management agents. Functions Select Event Management and System Management Monitoring). Click the 'GO' button to start the Sytem Management Installation/configuration.
The "System Management" installation / configuration takes ~25 minutes, please wait.	
Fig. 52: CA Check if installation was successful	The installation was successful, when the status has changed to partly running and the message "shutdown pending" is displayed in the in status bar of the Service UI.
System management installation / configuration pop-up will be displayed at the end of the installation.	
Fig. 53: CA Check if installation was successful	After the system re-boot: Check that: • Event Management has the status running • System Management Monitoring has the status enabled

System Management Functions

NOTE

If only one or two checkboxes are marked, CA Unicenter (System Management) is only partly installed. In this case, deinstallation of CA Unicenter (System Management) is required!

NOTE

If Auto Transfer is disabled and Event Management is selected, the collected data will not be transferred to the System Management unless a manual transfer is performed. This may result in data loss since the buffer size is limited (ring buffer).

Troubleshooting and useful information

System Management Installation failed

In case that the System Management Installation fails perform a manual de-installation as described below:

- Login as user administrator
- click **Start, Run** cmd
- **cd C:\MedCom\service\mwtools**
- **sysmgmt_install.exe -r**

- Check the status of CA Unicenter (System Management) by using the following command in the cmd tool window:

```
sysmgmt_install.exe -c<ENTER>
```

The following shows a successful deinstallation of CA UniCA Unicenter (System Management)

```
C:\MedCom\service\mwtools>sysmgmt_install.exe -c
```

```
*****
```

```
* Systemmanagement Managed Node Package - Version VE10D
```

```
*****
```

```
OS=WinXP/.Net
```

```
-----
```

```
Check Installation-Account
```

```
Account=administrator
```

```
-----
```

```
Check Registry-Keys
```

```
All Registry-Keys are present
```

```
-----
```

```
addr #0 is [192.168.1.1]
```

```
addr #1 is [192.168.2.1]
```

```
*****
```

```
* Systemmanagement Managed Node Package - Version VE10D Status:
```

```
*****
```

```
Backup Management: not installed
```

```
Event Management: not installed
```

```
Asset Agent not installed
```

```
Event Agent not installed
```

```
Hardware Agent not installed
```

```
Log Agent not installed
```

```
Software Delivery Agent not installed
```

```
Software Transport Agent not installed
```

```
System Agent not installed
```

```
Performance Management: not installed
```

```
Network Management: not installed
```

```
Rules/Filters Version: not installed
```

```
C:\MedCom\service\mwtools>
```

- If any entry from the "System Management Managed Node Package" shows installed, the program "`sysmgmt_install.exe -r`" must be run again. The status of CA Unicenter (System Management) has to be checked as well after the reboot.

If it still fails, contact the SRS help desk (Tel. xx49-9191-188080-191)

For problems within the US, call the US Uptime Service Center / technical support group.

For further investigation get the following logfile `C:\sysmgmt\EvtMgmtINIT.log` and `C:\TEMP\sysmgmt_install.log`.

Re installation of software fails

If the SW re installation fails e.g. due to a interrupt during the initial installation part of CD 1/2. The install software was partly portated to the system disk. At the next re installation, the install SW supposes that the SW installation is already performed with a dead lock as result.

Formatting the system disk

Work around:

- Insert the CD entitled "syngo MR CD 1/2" and reboot the system. Press "Alt-Ctrl-Del" simultaneously
- When the screen with the error message is displayed, press "CTRL-C" to cancel the installation.

Output	Input	Note
The DOS-Prompt appears on the left side of the screen.	a: ↵	If the drive letter is "C:\", then change to drive A:\ MS DOS including FDISK is loaded on RAM during the boot from CD 1/2. The ram content simulates a DOS floppy . No floppy is needed!
A:\	cd utils ↵	
List of options ...	fdisk ↵	Starts the partitioning program.

- Select the 4-th menu item, to get an overview of the existing partitions. Check now the file system type of the first partition.
- Press "**ESC**" to leave the current menu.
- Choose menu item "**3**" to delete Partitions.
- Depending on partition type, choose one of the items 1 to 4 to delete your first partition. In case you typed in a false menu item, "ESC" will bring you back to Step 8. This also applies, if message "NO PARTITION TO DELETE" appears after choosing a menu item.
- When you choose the right type for your first partition, the warning "DATA IN THE ... PARTITION WILL BE LOST" and "WHAT ...PARTITION DO YOU WANT TO DELETE [1]" appears
- Press "ENTER (↵)".

If a request for the partition name is asked, type in "MED_SYSTEM", which will also be displayed at the top of the screen.

- Continue with "Y" for "Yes"

NOTE

For some keyboard layouts the "Z" key will produce the "Y" in the DOS mode.

The message " ...PARTITION DELETED" appears.

- Exit the fdisk program by pressing "ESC" twice, press any key to reboot.

Try again to re install software, make sure that the "syngo MR CD 1/6" is still inserted in the CD-ROM drive.

System language and weight

Output	Input	Note
 Fig. 54:	Click Options Click Configura-tion Click Regional and Language Settings The "Regional Options" tab must be selected Select system language in the pull-down menu. Click Apply, Click OK	The system language must be set via the customer configuration platform within the MR application Patient weight unit for the "Patient Registration" input field, [lb] or [kg] is in relation to selected system language. English US = 1ball others = kg) Don't set the checkbox: "Set as system default locale" The language settings will be activated after following application restart
	Reboot the system	

Keyboard

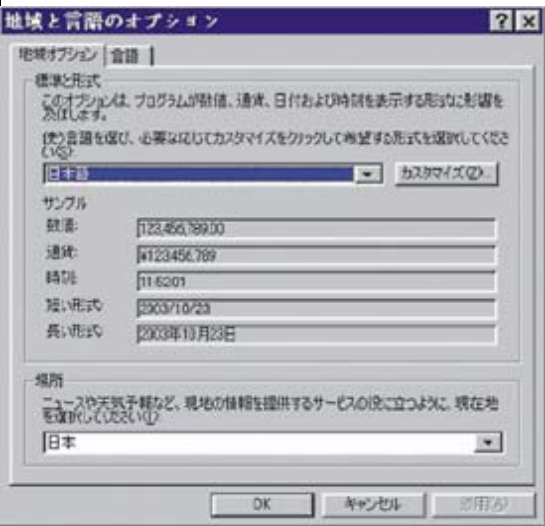
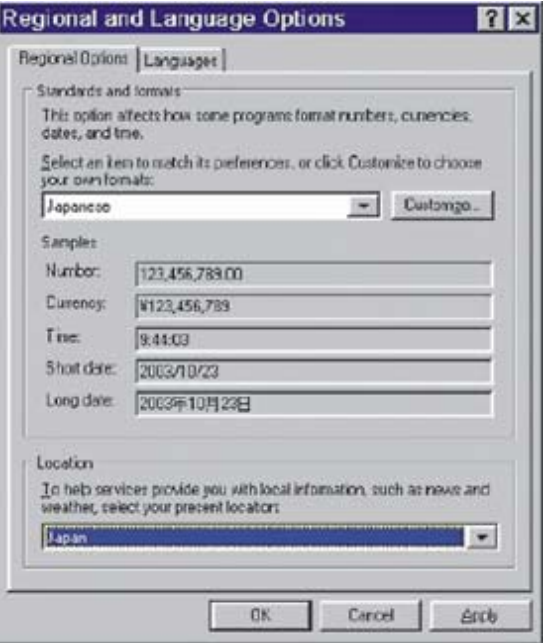
Refer to [\(Keyboard settings and Digit grouping symbol for meduser / p. 62\)](#).

Language Settings Japanese

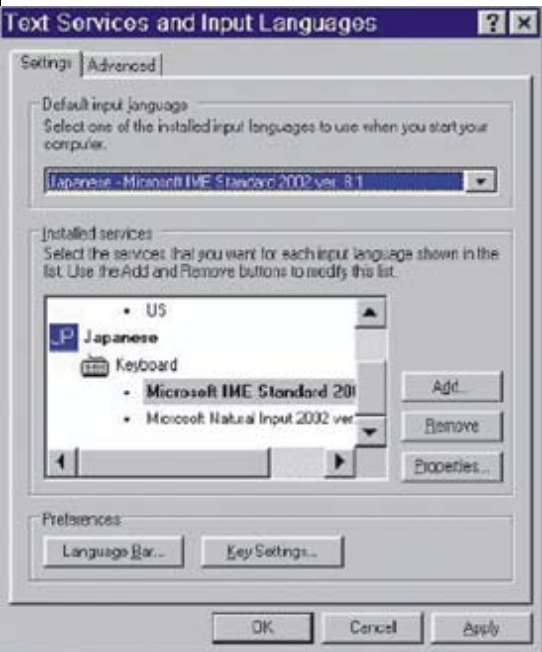
Switch language from international (english-US) to japanese

NOTE

The following steps can only be done, if the japanese version of WinXP (syngoMR 2004A Customer CD6/6) was used for the installation.

ACTION		MRC	MRSC
• Select Options -> Configuration in the main menu of the Meduser platform		☐	☐
• Double-click the Regional and Language Settings icon		☐	☐
• In the highlighted field choose Output popup in Japanese, if it was already switched during Windows XP Installation. 	Else the international view of the popup. 	☐	☐
Fig. 55: the first Item (see below)  (Hint: scroll up, if the following symbol is not present)	Fig. 56:		
Note: Be sure to use the long format for the date. It must look like yyyy/mm/dd or yyyy-mm-dd. If it is not set you have to change it. For it press the Customize button and select the correct setting for the date.			

ACTION		MRC	MRSC
Switch to the 2. tabcard “language” and choose the same Japanese symbol as in the first tabcard. Japanese symbol:  Fig. 57:			
Japanese view	or international (english-US) view	☐	☐
 Fig. 58:	 Fig. 59:		
Click on the "Details" Button and make following settings for your Keyboard:			

ACTION		MRC	MRSC
Japanese view	or international (english-US) view	☐	☐
 Fig. 60:	 Fig. 61:		
Click the Language Bar Button under preferences and select the two checkboxes like shown below:			
Japanese view	or international (english-US) view	☐	☐
 Fig. 62: Press the OK button.	 Fig. 63: Press the OK button.		
First press OK, next Apply. Note: Pressing "OK" in the "Regional and Language Options" Tab Card starts an automatic reboot, wait two steps with this action.		☐	☐
Open the Windows Explorer (It could be necessary to switch to "Advanced User" to open an Explorer)		☐	☐
Delete the file C:\MedCom\config\Mri\MRSystem\MrVCheck.ref		☐	☐
Press the OK Button and accept all Messages, which will be shown with "YES" or "OK". A Reboot follows!		☐	☐

ACTION	MRC	MRSC
Login as meduser (if autologin is not enabled) after the reboot and wait until the main UI is up.	☐	☐
IME settings: (by right mouse-click on "Input Mode"). 	☐	☐
<i>Fig. 64:</i> - Be sure that in the Configuration (Local service) the Dicom Character Set is set to "ISO 2002 IR13" for Japanese! - Start local service - Go in Configuration - Select under DICOM the link Character Set - Choose now the radio button "allow multiple character sets" - Select the Character set "ISO 2002 IR13" - Press "Save" and confirm the popup. - Restart the system by pressing "Home". Confirm the restart popups.	☐	☐

Switch language from japanese to international (english-US)

ACTION	MRC	MRSC
Select Options -> Configuration in the main menu of the Meduser platform	☐	☐
Double-click the Regional and Language Settings icon	☐	☐

ACTION	MRC	MRSC
Change on the first tabcard (Regional Options) the settings like in the picture below English - US.  Fig. 65:	☐	☐
Change on the second tabcard (Language) the settings like in the picture below English – US.  Fig. 66:	☐	☐
First press OK, next Apply. Hint: Pressing “OK” in the “Regional and Language Options” Tab Card starts an automatic reboot, wait two steps with this action.	☐	☐
Open the Windows Explorer (It could be necessary to switch to “Advanced User” to open an Explorer)	☐	☐

ACTION	MRC	MRSC
Delete the file C:\MedCom\config\Mri\MRSystem\MrVCheck.ref	☐	☐
Press the OK Button and accept all Messages, which will be shown with "YES" or "OK". A Reboot follows!	☐	☐

Language Settings Chinese

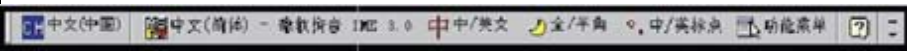
Switch language from international (english-US) to chinese

Note: The following steps can only be done, if the chinese version of WinXP (syngoMR 2004A Customer CD5/6) was used for the installation

ACTION		MRC	MRSC
Select Options -> Configuration in the main menu of the Meduser platform		☐	☐
Double-click the Regional and Language Settings icon		☐	☐
In the highlighted field choose Output popup in Chinese, if it was already switched during Windows XP Installation.		☐	☐
Else the international (english-US) view of the popup. Select the below shown entries.  Fig. 68:			
the following Item (see below)  Fig. 67:			
the following Item (see below)  Fig. 68:			
the following Item (see below)  Fig. 67:			
Switch to the 2. tabcard "language" and choose the same Chinese symbol as in the first tabcard. Chinese symbol:  Fig. 69:		☐	☐

ACTION		MRC	MRSC
Chinese view	or international (english-US) view		
 Fig. 70:	 Fig. 71:		
- Click on the "Details" Button Note: If the necessary Keyboard Layout/IME is not present, press the “Add” Button and select (see following popup) and press OK.  Fig. 72: - Make than the following settings for your Keyboard:	☐ ☐		

ACTION		MRC	MRSC
Chinese view  Fig. 73:	or international (english-US) view  Fig. 74:		
Click the Language Bar Button under preferences and select the two checkboxes like shown below:		☐	☐
Japanese view  Fig. 75: Press the OK button.	or international view  Fig. 76: Press the OK button.		
First press OK, next Apply. Note: Pressing "OK" in the "Regional and Language Options" Tab Card starts an automatic reboot, wait two steps with this action.		☐	☐
Open the Windows Explorer (It could be necessary to switch to "Advanced User" to open an Explorer)		☐	☐
Delete the file C:\MedCom\config\Mri\MRSystem\MRVCheck.ref		☐	☐
Press the OK Button and accept all Messages, which will be shown with "YES" or "OK". A Reboot follows!		☐	☐

ACTION	MRC	MRSC
Login as meduser (if autologin is not enabled) after the reboot and wait until the main UI is up.	☐	☐
You can switch the output language of IME to Chinese by a mouse-click on it: (See the red symbol in the middle).  Fig. 77: or to English  Fig. 78:	☐	☐

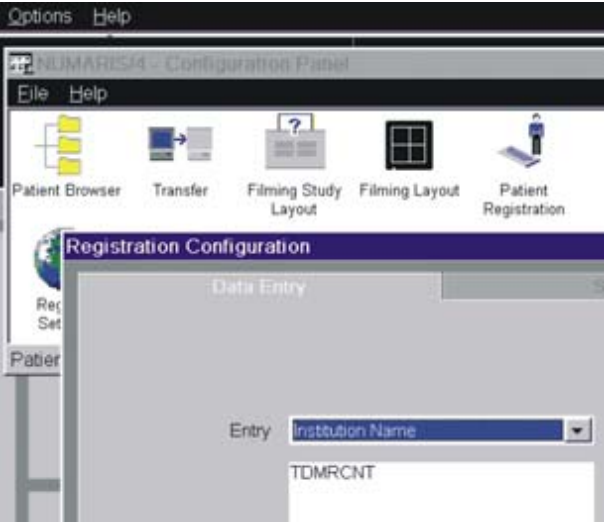
Switch language from chinese to international (english-US)

ACTION	MRC	MRSC
Select Options -> Configuration in the main menu of the Meduser platform	☐	☐
Double-click the Regional and Language Settings icon	☐	☐
Change on the first tabcard (Regional Options) the settings like in the picture below English - US.  Fig. 79:	☐	☐

ACTION	MRC	MRSC
Change on the second tabcard (Language) the settings like in the picture below English – US.  Fig. 80:	☐	☐
First press OK, next Apply. Hint: Pressing “OK” in the “Regional and Language Options” Tab Card starts an automatic reboot, wait two steps with this action.	☐	☐
Open the Windows Explorer (It could be necessary to switch to “Advanced User” to open an Explorer)	☐	☐
Delete the file C:\MedCom\config\Mri\MRSystem\MrVCheck.ref	☐	☐
Press the OK Button and accept all Messages, which will be shown with “YES” or “OK”. A Reboot follows!	☐	☐

Institution name

The Customer UI must be up and running.

Output	Input	Note
 Fig. 81:	Click Options Click Configuration Click Patient Registration Click Institution Name	The Institution name will be part of the DICOM header

Regional settings

Date and Time syntax

To configure the date and time syntax, proceed as follows: select Options, Configuration, Regional Settings, **Time or Date** in the MR application. It is recommended to select a 4 digit syntax for the year in the "Short date style" pull-down menu to ensure the correct patient age on film printouts (e.g., for patients born in the 19th century or infants).

When you are finished defining the regional settings, restart the MR system to activate the changes.

System Manager

Restart Application

In case the MR application has crashed or is frozen, try starting the System Manager as follows:

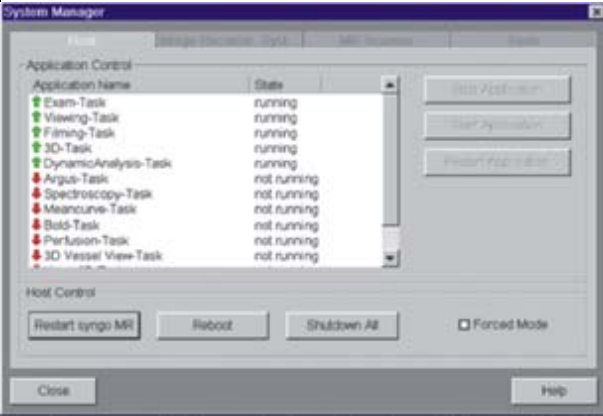
Output	Input	Note
	Simultaneously press [Ctrl]-[Alt]-s	
	Click Restart syngo MR	

Fig. 82:

Copy SIEMENS Programs to the customer tree

Use the following procedure to copy SIEMENS applications to the customer tree. The customer can modify these templates according to his own requirements.

- Start the Service platform

Output	Input	Note
	cmd.exe	
	InstallProtocols	

NOTE

If there are existing protocols in the customer tree, the "InstallProtocols" command creates a new tree without overwriting the existing one. The new SIEMENS program folders are numbered "001", "002", etc.

IP Security

The syngo-MR systems are protected with a service "ipsec".

Other network nodes, like DICOM nodes, can only access the DICOM services that your console provides if they are included in the "Generic Transport Filter" of ipsec.

To list nodes that are allowed to access your console, click **Reports**; select **IP Security** below **SiteInfo**.

Adding host to IP security list

To add a node manually into the IP security list, use following procedure:

At the "MRC", simultaneously press **Ctrl - Esc**. Select **Start, Advanced User**; enter the advanced user password; click **OK**

Simultaneously press **Ctrl - Esc**, click **Start, Run cmd** enter the command: **runas /user:administrator cmd** enter the password

In the newly started cmd tool you have administrator rights.

(replace the address "146.254.222.333" in the following **example** with the IP address of your remote node)

Input in the cmd window	Output
ipseccmd -f (0+146.254.222.333)	The command completed successfully

Deleting all IP security settings

To delete all nodes manually from the IP security list, use following procedure:

At the "MRC", simultaneously press **Ctrl - Esc**. Select **Start, Advanced User**; enter the advanced user password; click **OK**

Simultaneously press **Ctrl - Esc**, click **Start, Run cmd** enter the command: **runas /user:administrator cmd** enter the password

Input in the cmd window	Output
ipseccmd -u	The command completed successfully

Adding all configured hosts to IP security list

To add configured nodes manually into the IP security list, use following procedure:

At the "MRC", simultaneously press **Ctrl - Esc**. Select **Start, Advanced User**; enter the advanced user password; click **OK**

Simultaneously press **Ctrl - Esc**, click **Start, Run cmd** enter the command: **runas /user:administrator cmd** enter the password

Input in the cmd window	Output
C:\MedCom\utils\SetIPSec.bat /install	... Setting dynamic security policy for details see C:\MedCom\log\SetIPSec.log

Context Vision

The Context Vision filter option requires an additional license key that is defined by the **dongle** and the system **serial number**.

- Start the service platform (service level 7 is required)

Requested data for cvlicense generation

Output	Input	Note
	Click Utilities Select Escape to OS	
 <i>Fig. 83: Get Host ID</i>	Type $\chi\omega\lambda\iota\chi\epsilon\nu\sigma\epsilon$ ↵	

Your notes:

Product	
Serial n/o	
Host Name:	
Dongle ID	Without 7- , 9- or Dallas dongles without 8-

All licenses are handled by TD HSC 1 in Germany (e-mail: TDHSC@med.siemens.de, phone: +49-9191-188080). (Please ask for the "call center".)

Activate license for Context Vision

Output	Input	Note
	Click Utilities Select Escape to OS	
 <i>Fig. 84: Enter license for Context Vision</i> Enter license for Context Vision	Type cvlicense [license key] ↵	

What to do next?		
If you are configuring the MRC or MRSC, refer to (Restoring site specific data MR 2004A Rev. B on MR 2004A Rev. B / p. 84) 	If you are installing / adding an MRSC, refer to (Import configuration data from CD-R to MRSC / p. 60) 	If you are performing the first configuration, refer to (Configuration / p. 88) 

MR auto login

Access NT login

Press SHIFT during the "blue screen" boot phase to access the NT login.

NOTE

To disable/enable MR autologin, refer to On-line Help texts "Configuration", "Users".

syngo startup / shutdown counter

How to terminate the syngo counter

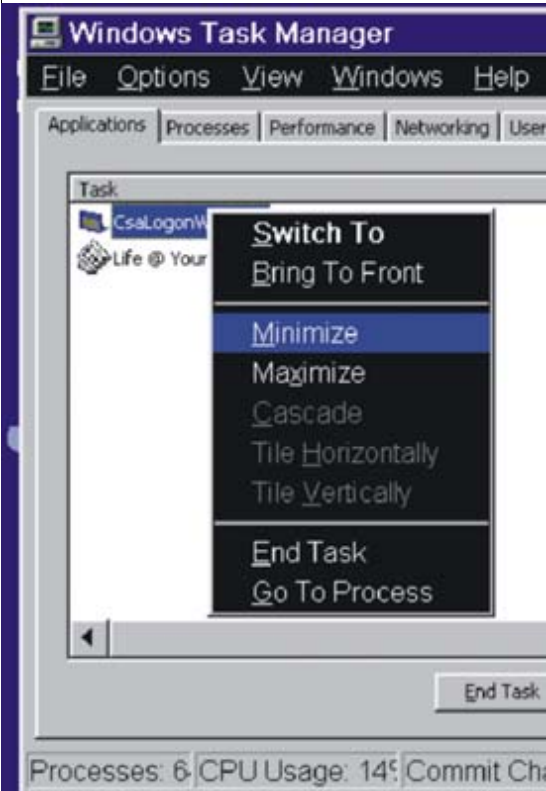
Output	Input	Note
 The screenshot shows the Windows Task Manager window with the 'Applications' tab selected. A right-click context menu is open for the 'CsaLogonWindow' task. The menu options are: Switch To, Bring To Front, Minimize (highlighted), Maximize, Cascade, Tile Horizontally, Tile Vertically, End Task, and Go To Process. The status bar at the bottom shows 'Processes: 6', 'CPU Usage: 14%', and 'Commit Ch'.	Terminating the syngo counter: Press Simultaneously Ctrl, Shift, and ESC Select the tab Applications right click CsaLogonWindow select Minimize	To view actions during startup and shutdown, you can minimize the syngo counter

Fig. 85: syngo counter

Terminate deadlock of crashed Service UI

Internet Explorer can cause a system deadlock during a service session.

How to terminate the Service IU (Internet Explorer)

Output	Input	Note
 Fig. 86: Windows NT Task Manager	Simultaneously press "Ctrl-Alt-Delete" Click the "Applications" tab Select the Siemens Med Service Software task in the task list Click "End Task"	The service platform crashes, the system screen is white (The User Interface will recover within one minute)

Expert - i

Following is only valid for systems with the "**Expert-i**" option To check if the "MR_Expert_i" licence is available, use the Service: UI, Configuration, Service, Licensing.

- Log in as **administrator** at the **MRC or MRSC**

Output	Input	Note
	Insert the "installtion CD" for the option "Expert-i" to the CD or DVD ROM-drive	
	Open the "Windows file Explorer"	
	Navigate to CD or DVD drive (R:)	
	Double-click the Host.exe	This installs the Expert-i option
	Double-click the Host.exe	Due to the MR firewall the skript must be run once more
	Re-boot, login as meduser	The Expert-i icon will be found in the status-bar

Measurement system configuration matrix

SW Product No.:	Product Name
024	Harmony
087L	Harmony Impact (Large cabine)
087S	Harmony Impact (Small cabine)
047	Symphony
089	Symphony Vision
085	Sonata
086	Sonata Vision
007	Concerto
083	Allegra
091	Trio

SW Product No.:	024	087L	087S	047	089	085	086	007	083	091
MPCU										
PPC/266 MHz	x	x	x	x	x	x	x	x	-	-
PPC/300 Mhz	x	x	x	x	x	x	x	x	x	x
RX channels										
RX4	x	x	x	x	x	x	x	x	x	-
2. RX4	o	o	o	o	o	o	o		-	x
Receiver										
Receiver 1	x	x	x	x	x	x	x	x	x	-
Receiver 2	o	x	x	o	x	x	x	o	o	-
Receiver 3	o	o	o	o	o	o	o		-	-
Receiver 4	o	o	o	o	o	o	o		-	x
Transmitter										
TRA1	x	x	x	x	x	x	x	x	x	x
TRA2				o	o	o	o	-	-	o
EFI										
No EFI_DAS	-	-	-	-	-	-	-	x	-	-
EFI_DAS	-	-	-	-	-	-		x	-	-
EFI_DAS with Sensor	-	-	-	-	-	-	-	x	-	-
Mains frequency										
50 Hz	-	-	-	-	-	-	-	x	-	-
60 Hz	-	-	-	-	-	-	-	x	-	-
Railway frequency (only valid for sites located close to a railway track)										
16 2/3 Hz	-	-	-	-	-	-	-	x	-	-
25 Hz	-	-	-	-	-	-	-	x	-	-
Gradient amplifier										
Platform/400V	x			x						
Platform/800V	x			x						
Cascade/500V	x			x						
Cascade/1000V 300A	x			x						
Cascade/1000V 380A	x			x						
Cascade/1500V 380A	x			x						

SW Product No.:	024	087L	087S	047	089	085	086	007	083	091
Cascade/2000V 380A	x	x	x	x	x					
Cascade/2000V 500A						x	x		x	x
GPA007								x		
Gradient coils										
AS39S	x			x						
AS39SR	x			x						
AS39T	x	x	x	x	x					
AS Sonata						x				
AS Sonata Vision							x			x
GC018								x		
AC 44	-							x	x	-
Stimulation monitor										
GPA_INT	x	x	x	x	x	x	x	-	x	x
dB/dt Look Ahead Guideline										
SAFE model	x	x	x	x	x	x	x		x	x
FDA	x	x	x	x	x	x	x		x	x
IEC	x	x	x	x	x	x	x		x	x
MHW	x	x	x	x	x	x	x		x	x
RF Amplifier										
DORA (RFPK2204)	x	x	x	x	x	x	x			
CORA (RFPK2137)				o	o	o	o	x		
RFPK (LPPA 13040)				o	o	o	o	-	x	-
RFPK (LPPA 13080)				-	-	-	-	-	-	x
RF Resonator										
BC LP rx/tx	x	x	x						-	-
BC CP rx/tx				x	x	x	x		-	-
BC 2 CP rx/tx				x	x	x	x		-	-
BC LP -/tx								x	-	-
RF safety watchdog										
RFSWD	x	x	x	x	x	x	x		x	x
Number of Body Tune Boxes										

SW Product No.:	024	087L	087S	047	089	085	086	007	083	091
One BTB installed	x	x	x						-	-
Two BTBs installed				x	x	x	x		-	-
BTB Types										
BTB 024	x								-	-
BTB 024S			x						-	-
BTB 024L		x							-	-
BTB 047 AS<5				x					-	-
BTB 047 AS>=5				x					-	-
BTB 047 Syngo				x	x				x	-
BTB Sonata						x			-	-
BTB Sonata Vision							x		-	-
BTB 063								x	-	-
BTB 047-2				x	x	x	x		-	-
BTB 3TWB									-	x
RFCI										
CAPTICON								x		
2.DYSCON	x	x	x	x	x	x	x		x	x
Receiver										
LCCS	x	x	x	x	x	x	x		x	x
COSIMUX								x		
Magnet										
OR41		x	x							
OR42		x	x							
OR35					x		x			
OR24	x									
OR70				x		x				
OR64										x
OR92									x	
Sumitomo								x		
Magnet Supervision										
MSU	x	x	x	x	x	x	x		x	x

SW Product No.:	024	087L	087S	047	089	085	086	007	083	091
TCU								x		
Active Shim										
Additional shim	o	o	o	o	o	x	x		x	x
array shim / standard	x	x		x	x	x	x		-	-
array shim / small cabin			x						-	-
mechanic shim								x	-	-
Front end	o	x	x	o	x	x	x	o	-	-
Patient table										
Distance Standard	x			x		x			-	
Distance Impact		x	x						-	
Distance Vision					x		x		-	x
Open-Table								x		
OP Table				o				o		
OHF Table										
LCD-Display+Banane	x	x	x	x	x	x	x			x
DTP								x		
Spectroscopy (structure related to appl. packages)										
Proton spectroscopy				o	o	o	o		x	x
SAMI-Ansteuerung				o	o	o	o			
2.RFPA				o	o	o	o			
Cooling system										
RCA 1 cooler	x	x	x	x	x		x		-	-
RCA 2 cooler						x			x	x
GRACCO								x		
Mobile										
Mobile enforce	o			o		o			-	-
Coils										
CP Flex large	o	o	o	o	o	o	o			
CP Flex small	o	o	o	o	o	o	o			
Flex loop large	o	o	o	o	o	o	o			
Flex loop small	o	o	o	o	o	o	o			

SW Product No.:	024	087L	087S	047	089	085	086	007	083	091
Loop medium										
Endorectal	o	o	o	o	o	o	o			
TMJ right	o	o	o	o	o	o	o			
TMJ left	o	o	o	o	o	o	o			
CP peripheral angio	o	o	o	o	o	o	o			
CP breast array	o	o	o	o	o	o	o			
CP wrist										
Shoulder array	o	o	o	o	o	o	o			
CP extremity	o	o	o	o	o	o	o			
CP body array	o	o	o	o	o	o	o			
CP neck array	o	o	o	o	o	o	o			
CP spine array	o	o	o	o	o	o	o			
Head array	o	o	o	o	o	o	o			
Head CP (TX/RX)	-	-	-	-	-	-	-	-	x	x
Cardiac array	o	o	o	o	o	o	o			
MP 11								o		
MP 16								o		
MP 21s								o		
MP 28								o		
Head neck array								o		
Body spine array								o		

o = Option

Dicom

Dicom and UNIX nodes

For UNIX-based SIEMENS systems, such as NUMARIS 3x nodes, MagicView 1000 etc., make sure to define separate logical names for each service, e.g. Storage or Query & Retrieve.

Dicom configuration data for SIEMENS products

	Product	DICOM Service	AET	port no.
CT	Somaris 4:	Storage SCU (send)	<PLA prefix> + DFOS	-
	Somatom Plus4	Storage SCP (receive)	<PLA prefix> + DIC1	50082
	Somatom AR (Spiral)	Query/Retrieve SCU	<PLA prefix> + DQRU	-
	Somaris 4.5:	Query/Retrieve SCP	<PLA prefix> + DQRY	50089
	Balance	Basic Print SCU	<PLA prefix> + DBPR	-
		Modality Worklist SCU	<PLA prefix> + DBWM	-
CT	Somaris 5:	Storage SCU (send)	STU_<hostname> ¹	-
	Somatom VZ	Storage SCP (receive)	STU_<hostname>	104
		Query/Retrieve SCU	STU_<hostname>	-
		Query/Retrieve SCP	STU_<hostname>	104
		Basic Print SCU	MERGE_PRINT_SCU	-
		Modality Worklist SCU	HRI_<hostname>	-
CT	3D Virtuoso VA30	Storage SCU (send)	<HOSTNAME> ²	-
		Storage SCP (receive)	<HOSTNAME>	4444
		Query/Retrieve SCU	<HOSTNAME>	-
		Query/Retrieve SCP	<HOSTNAME>	4444
		Basic Print		

1. These are the default values for the AET. It is recommended to use the <hostname> for the AET. Future software versions will change the default value to the <hostname> as well.
2. Host name in uppercase letters

	Product	DICOM Service	AET	port no.
MR	Numaris 3: Magnetom Impact Magnetom Expert Magnetom Vision Magnetom Open	Storage SCU (send)	<PLA prefix> + DFOS	-
		Storage SCP (receive)	<PLA prefix> + DIC1	50082
		Query/Retrieve SCP	<PLA prefix> + DQRY	50089
		Basic Print SCU	<PLA prefix> + DBPR	-
		Modality Worklist SCU	<PLA prefix> + DBWM	-
MR	Numaris 3.5: Harmony Symphony			
MR	syngo MR (Numaris 4)	Storage SCU (send)	STU_<hostname> ¹	-
		Storage SCP (receive)	STU_<hostname>	104
		Query/Retrieve SCU	STU_<hostname>	-
		Query/Retrieve SCP	STU_<hostname>	104
		Basic Print SCU	MERGE_PRINT_SCU	-
		Modality Worklist SCU	HRI_<hostname>	-

1. These are the default values for the AET. It is recommended to use the <hostname> for the AET. Future software versions will change the default value to the <hostname> as well.

	Product	DICOM Service	AET	port no.
SHS	MagicView 1000 VA31 / VB32	Storage SCU (send)	<PLA prefix> + DRSP ¹	-
		Storage SCP (receive)	<PLA prefix> + DRSP	50082
		Query/Retrieve SCU	<PLA prefix> + DQRY	-
		Query/Retrieve SCP	<PLA prefix> + DQRY	50089
		Basic Print SCU	<PLA prefix> + DBPR	-
SHS	MagicStore	Storage SCU (send)	<PLA prefix> + DRSP	-
		Storage SCP (receive)	<PLA prefix> + DRSP	50082
		Query/Retrieve SCP	<PLA prefix> + DQRY	50089 ²
		Query/Retrieve SCP	<PLA prefix> + DPQR	50090 ³
		Basic Print SCU	<PLA prefix> + DBPR	-

	Product	DICOM Service	AET	port no.
SHS	MagicView 50 / 300	Storage SCU (send)	configurable ⁴	-
		Storage SCP (receive)	configurable	104
		Query/Retrieve SCU	configurable	-
SHS	MagicLink D / MagicServe	Storage SCU (send)	<PLA prefix> + DFOS	-
		Storage SCP (receive)	<PLA prefix> + DIC1-8 ⁵	50082

1. Configurable with login: install -edit -config
2. Short Term database (IMS or IM)
3. Long Term database (PDIR or MM)
4. it is recommended to use the <hostname> as AET
5. 8 import channels (DIC1 to DIC8) with different PACSnet destinations can be configured

Migration from Frame-Word to CMS.

Chapter 5: Information added how to Backup/Restore site specific data on USB device.

There are no Hazard IDs in this document.

