

SIGNA® HORIZON™ 1.5T & 1.0T UFI–COMPATIBLE SURFACE COILS WARNING LABEL INSTRUCTIONS

FOR THE GE SERVICE ENGINEER:

Product Description:

The following procedures need to be followed for all new and upgrade 1.5T Signa Horizon system installations.

Note

The Surface Coil UFI–compatibility issue and warning labels plan do **NOT** apply to 1.0T Signa Horizon systems.

With the purchase or upgrade of a 1.5T Horizon, HiSpeed, or EchoSpeed, there are several RF coils that will no longer be compatible with the Signa Horizon system. It is required that these non–compatible coils be removed from the site and **returned** to GE. A new coil, of the same type, will be supplied for each of these non–compatible coils before they are removed from the site. This will be accomplished by a marketing promotion, special no–cost S–CAT's, or FMI.

Note

These procedures **do not** apply to 1.5T Signa “Advantage” 5.5 Software Upgrade only systems. For these systems, selecting the “Advantage” Gradient Configuration file option during Load from Cold, will set B1 RMS to the appropriate value.


For NON–GE Approved surface coils: Refer the customer to the original coil manufacturer for information on coil compatibility with Signa Horizon.

Note

GE makes no claims as to compatibility of coils not approved by GE for use on Signa systems. We have no way of knowing if any of these non–GE coils are compatible with 5.5 software and/or hardware. The customer uses these coils at their own risk, and GE accepts no responsibility for their use.

Please return unused coil labels to Robert Jacobs, W–827. For more information, please contact me at (414) 521–6733 (Email JACOBS–R).

Service Readiness:

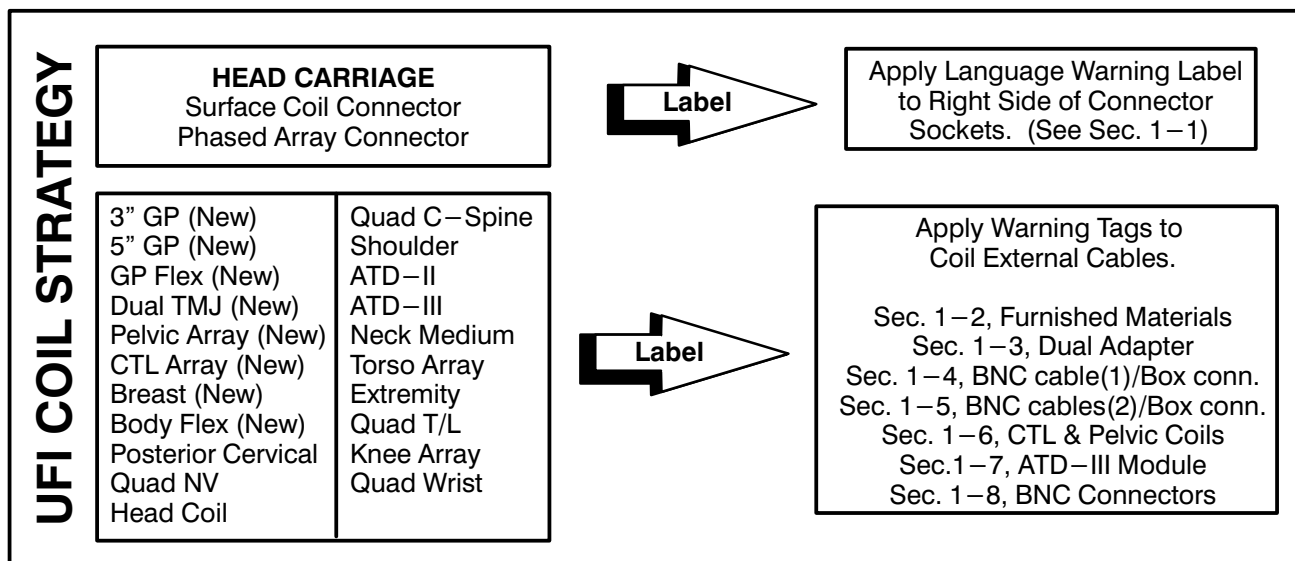

(Continued)



SIGNA® HORIZON™ 1.5T & 1.0T UFI-COMPATIBLE SURFACE COILS WARNING LABEL INSTRUCTIONS

1-1 ADD WARNING LABELS TO HEAD CARRIAGE TO INDICATE UFI-COMPATIBILITY

It is necessary to add special warning labels to the surface coil and phased array connectors on the Signa Horizon 1.5T Head Carriage. There are eight different language warning labels provided, English, French, Spanish, German, Italian, Portuguese, Japanese, and Chinese. (In most cases, the English language labels will be on the delivered Head Carriage.) Add the appropriate language warning labels as shown:

Labeling Instructions For Head Carriages:

1. Clean surface with a damp cloth to remove dirt and oil. Allow time to dry.
2. Remove label backing. Attach label to the right of head & phased array connectors at the head carriage. Apply foreign language labels on top of existing English language labels, if necessary.
3. Allow labels to cure for 12 hours with minimum handling to maximize adhesive strength.

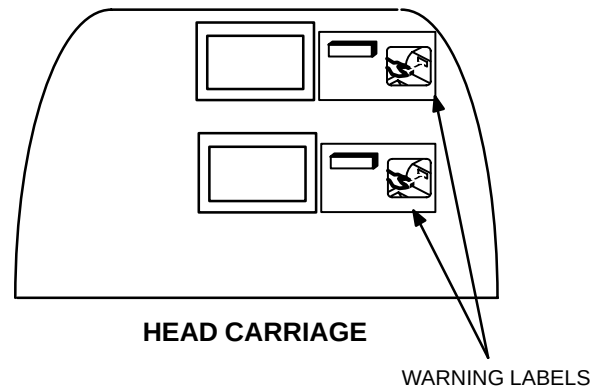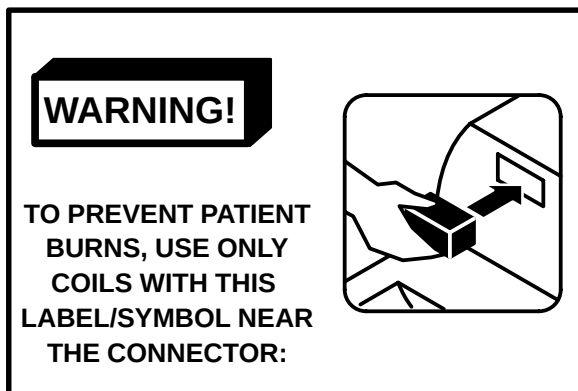

1-2 FURNISHED MATERIALS

A kit, 2146243 (FRU-2), is attached to inside front frame of the RF/PEN Cabinet above the RF System Controller Module. This kit includes:

Item	Part Number	FRU	Name	Quantity	Description (Remarks)
1	2145936	X	INSTALLATION	1	INSTALLATION DIRECTIONS
2	2146833	2	VINYL STRIPS	20	WHITE VINYL STRIP, 1" X 5", ADHESIVE BACK
3	2146832	2	WARNING TAG	20	UFI WARNING TAG, ADHESIVE BACK
4	2146831 THRU -8	2	WARNING LABELS	2 EA.	UFI WARNING LABEL, ADHESIVE BACK, 8 LANGUAGES

Service Features:

1-3 LABELING INSTRUCTIONS FOR COILS WITH DUAL ADAPTER

The Following Coils Are UFI-compatible. These coils should be labeled based on the coils attached:

COIL (1.5T)	UFI M-CAT#	UFI GE #	LABEL ATTACHMENT
1.5T Dual Adapter	M1087DA	N/A	Refer to Coil Attached.

(Continued)

SIGNA® HORIZON™ 1.5T & 1.0T UFI-COMPATIBLE SURFACE COILS WARNING LABEL INSTRUCTIONS

1-4 LABELING INSTRUCTIONS FOR COILS WITH ONE BNC CABLE/BOX CONNECTOR



The Following Coils Are UFI-compatible. These coils need to have UFI-compatible warning tags attached:

COIL (1.5T)	UFI M-CAT#	UFI GE #	LABEL ATTACHMENT
1.5T ATD-II Coil	M1085EE	46-317599P1	Label Tag (Corrugated tube)
1.5T Breast Coil	M1085BR (After FMI 60413)	46-320785P2	Label Tag (Corrugated tube) See Note 1.
1.5T Knee Array Coil	M0187AA	46-320406P1	Label Tag (Corrugated tube)
1.5T Quad NV	M1085HN	2118138	Label Tag (Corrugated tube)
1.5T Quad Wrist	E8801K	MAI # 10155	Label Tag (Corrugated tube)
NOTE 1: Apply Label to phased array (46-320756P1) external cable.			

Labeling Instructions For BNC Cable with Box Connector:

- Clean surface of cable with a damp cloth to remove dirt and oil. Allow time to dry.
- Remove backing of white vinyl strip. Position vinyl strip as close to connector as possible without interfering with connector, strain relief, and existing cable labeling. (See Illustration below right.)
- Wrap entire length of white vinyl strip as tight as possible around cable jacket.
- Attach Warning Tag as follows: (See Illustration below left.)
 - Remove backing of clear strip of the warning tag, apply clear strip to vinyl. (Keep warning label parallel to connector.)
 - Wrap around white vinyl strip one time until strip is facing warning label.
 - Remove backing of warning label and apply clear strip to surface of warning label.
 - Wrap clear strip around edge of warning label.
 - Apply clear strip to other surface of warning label.
 - Cut excess clear strip where clear strip meets white vinyl strip.
- Allow label to cure for 12 hours with minimum handling to maximize adhesive strength.

**NOTE LABEL ORIENTATION
PARALLEL TO HOUSING.**

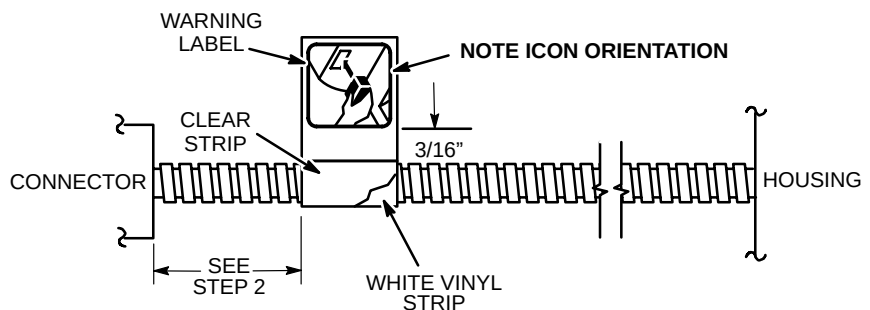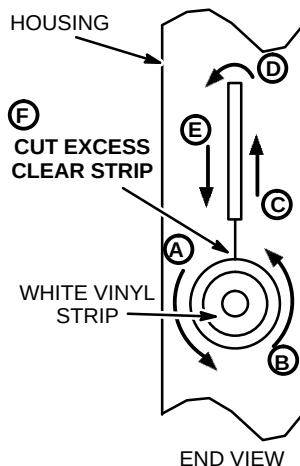

ATD-II, BREAST (PA), KNEE ARRAY, QUAD NV, QUAD WRIST

(Continued)

SIGNA[®] HORIZON[™] 1.5T & 1.0T UFI-COMPATIBLE SURFACE COILS WARNING LABEL INSTRUCTIONS


1-5 LABELING INSTRUCTIONS FOR COILS WITH TWO BNC CABLES/BOX CONNECTOR

The Following Coils Are UFI-compatible. These coils need UFI-compatible warning tags attached:

COIL (1.5T)	UFI M-CAT#	UFI GE #	LABEL ATTACH.
1.5T Torso Array Coil	M1087TA	2104700	Label Tag (Corrugated tube)
1.5T Breast Coil	M1085BR (After FMI 60413)	46-320785P2	Label Tag (Corrugated tube) See Note 1.
1.5T Head Coil	46-282118G2	46-282118G2	Label Tag (BNC Cable)
NOTE 1: Apply Label to bilateral (46-320765P1) external cables.			

Labeling Instructions For Two BNC Cables with Box Connector:

- Clean surface of cable with a damp cloth to remove dirt and oil. Allow time to dry.
- Remove backing of white vinyl strip. Position vinyl strip as close to connector as possible without interfering with connector, strain relief, and existing cable labeling. (See Illustration below.)
- Wrap entire length of white vinyl strip as tight as possible around cable jacket.
- Attach Warning Tag as follows: (See Illustration below left.)
 - Remove backing of clear strip of the warning tag, apply clear strip to vinyl. (Keep warning label parallel to connector.)
 - Wrap around white vinyl strip one time until strip is facing warning label.
 - Remove backing of warning label and apply clear strip to surface of warning label.
 - Wrap clear strip around edge of warning label.
 - Apply clear strip to other surface of warning label.
 - Cut excess clear strip where clear strip meets white vinyl strip.
- Repeat for Head Coil Only:** Repeat the above process to add second label tag as close to housing as possible without interfering with strain relief and existing cable labeling. **Reverse direction of wrapping in step 4.** (Keep warning label parallel to housing.)
- Allow label to cure for 12 hours with minimum handling to maximize adhesive strength.

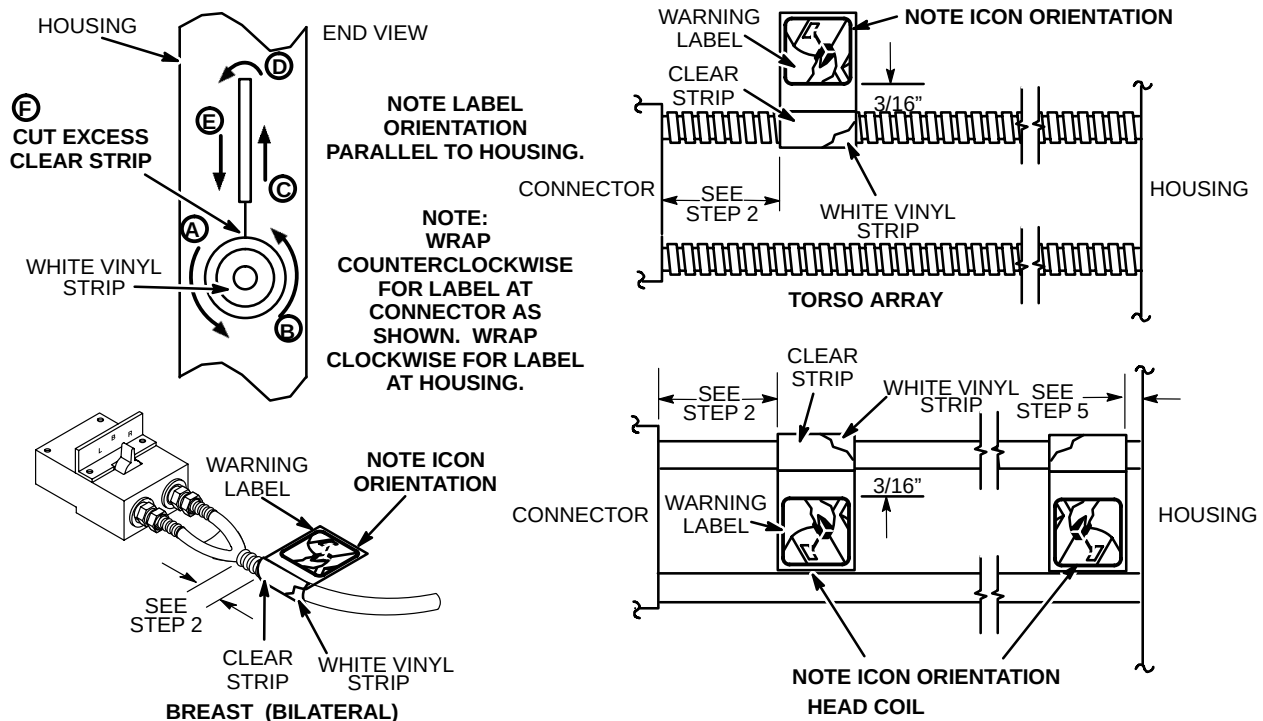

(Continued)

SIGNA® HORIZON™ 1.5T & 1.0T UFI-COMPATIBLE SURFACE COILS WARNING LABEL INSTRUCTIONS


1-6 LABELING INSTRUCTIONS FOR CTL & PELVIC COILS

The Following Coils Are UFI-compatible. These coils need to have UFI-compatible labels attached: (These coils use flat external cables.) These coils are being manufactured with the warning tag applied, however, please verify.


COIL (1.5T)	UFI M-CAT#	UFI GE #	LABEL ATTACHMENT
1.5T CTL Array Coil	M1087CT	2132617	Label Tag (Flat cable)
		2104799	Label Tag (Flat cable)
1.5T Pelvic Array Coil	M1087PA	2144100	Label Tag (Flat cable)

Labeling Instructions For CTL and Pelvic Array Coils:

1. Clean surface of flat cable with a damp cloth to remove dirt and oil. Allow time to dry.
2. Remove warning label & clear strip backings. Attach warning label to top of flat cable as shown:
 - A. Wrap clear strip around back of flat cable,
 - B. Wrap clear strip over label,
 - C. Cut excess clear strip on back side of flat cable.
4. Allow label to cure for 12 hours with minimum handling to maximize adhesive strength.

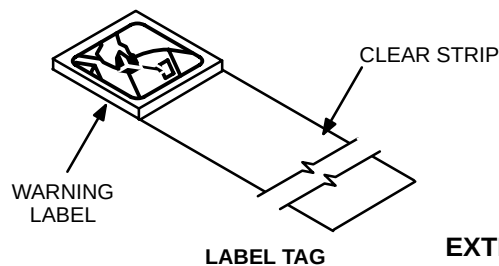

EXTERNAL FLAT CABLES
CTL, PELVIC COILS

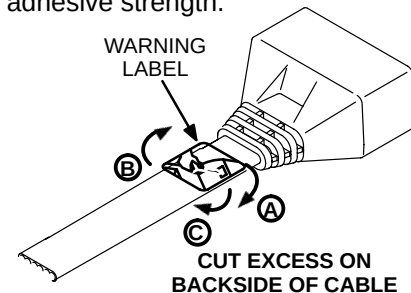


1-7 LABELING INSTRUCTIONS FOR ATD-III MODULE

The Following Coil Is UFI-compatible. This coil needs to have a UFI-compatible label attached:


COIL (1.5T)	UFI M-CAT#	UFI GE #	LABEL ATTACHMENT
1.5T ATD-III Coil	M1087AT	46-317600P1	Label Tag (ATD-III Housing)

Labeling Instructions For ATD-III Module:

1. Clean surface of module with a damp cloth to remove dirt and oil. Allow time to dry.
2. Cut clear strip from warning label as shown below:
3. Remove warning label backing. Attach label to top of ATD-III Module as shown:
4. Allow label to cure for 12 hours with minimum handling to maximize adhesive strength.

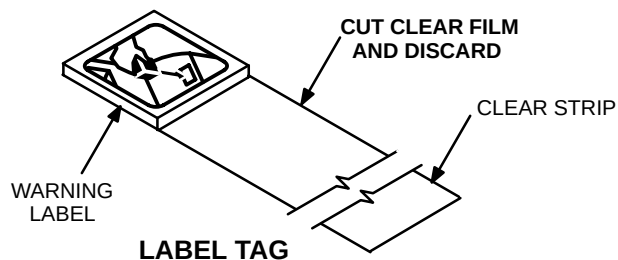

ATD-III MODULE

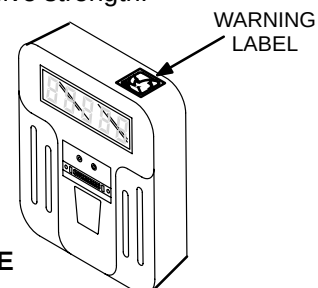

(Continued)

SIGNA® HORIZON™ 1.5T & 1.0T UFI-COMPATIBLE SURFACE COILS WARNING LABEL INSTRUCTIONS


1-8 LABELING INSTRUCTIONS FOR COILS WITH BNC CABLE/BNC CONNECTOR

The Following Coils Are UFI-compatible. These coils need to have a warning label at the external cable connector **and** an additional warning label attached on the external cable near the coil housing: (This is to address usage of these coils with Dual Array Adapter.)

COIL (1.5T)	UFI M-CAT#	UFI GE #	LABEL ATTACHMENT
1.5T GP Flex Coil	M1085GR	2128554	Label Tag (Corrugated tube)
1.5T Neck-Medium	M1085AK	46-287396P10	Label Tag (Corrugated tube)
1.5T Posterior Cervical	M1085AM	46-287807P12	Label Tag (Corrugated tube)
1.5T Q-T/L-Spine w/P	M1085AP	46-307387P10	Label Tag (Corrugated tube)
1.5T Extremity Surf	M1085AC	46-265442G1	Label Tag (Corrugated tube)
1.5T Quad C-Spine	M1085AR	46-320422P1	Label Tag (Corrugated tube)
1.5T Shoulder Coil	M1085AN	46-301058P9	Label Tag (Corrugated tube)
1.5T 3" Round GP	M1085GA	2127315	Label Tag (Corrugated tube)
1.5T 5" Round GP	M1085HA	2127316	Label Tag (Corrugated tube)
1.5T Body Flex	M1085BE	2148690	Label Tag (Corrugated tube) (On both coil external cables.)

Labeling Instructions For BNC Cable with BNC Connector:

- Clean surface of cable with a damp cloth to remove dirt and oil. Allow time to dry.
- Remove backing of white vinyl strip. Position vinyl strip as close to connector as possible without interfering with connector, strain relief, and existing cable labeling. (See Illustration below.)
- Wrap entire length of white vinyl strip as tight as possible around cable jacket.
- Attach Warning Tag as follows: (See Illustration below left, **wrap counterclockwise at connector.**)
 - Remove backing of clear strip of the warning tag, apply clear strip to vinyl. (Keep warning label parallel to connector.)
 - Wrap around white vinyl strip one time until strip is facing warning label.
 - Remove backing of warning label and apply clear strip to surface of warning label.
 - Wrap clear strip around edge of warning label.
 - Apply clear strip to other surface of warning label.
 - Cut excess clear strip where clear strip meets white vinyl strip.
- Repeat the above process to add second label tag as close to housing as possible without interfering with strain relief and existing cable labeling. (See Illustration below right, **wrap clockwise direction in step 4 at housing.** (Keep warning label parallel to housing.)
- Allow label to cure for 12 hours with minimum handling to maximize adhesive strength.

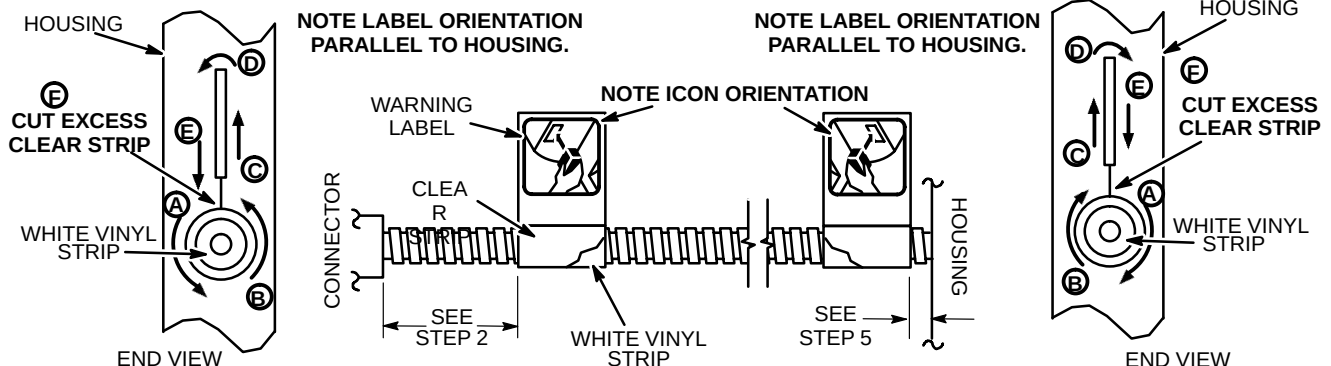

GP FLEX, NECK-MEDIUM, POSTERIOR CERVICAL, Q-T/L-SPINE, EXTREMITY, QUAD C-SPINE, SHOULDER, 3" GP, 5" GP, BODY FLEX