

Product Description:

The following procedures need to be followed for all new and upgrade 1.5T Signa Horizon system installations.

Note

The Surface Coil UFI-compatibility issue and warning labels plan do **NOT** apply to 1.0T Signa Horizon systems.

With the purchase or upgrade of a 1.5T Horizon, HiSpeed, or EchoSpeed, there are several RF coils that will no longer be compatible with the Signa Horizon system. It is required that these non-compatible coils be removed from the site and **returned** to GE. A new coil, of the same type, will be supplied for each of these non-compatible coils before they are removed from the site. This will be accomplished by a marketing promotion, special no-cost S-CAT's, or FMI.

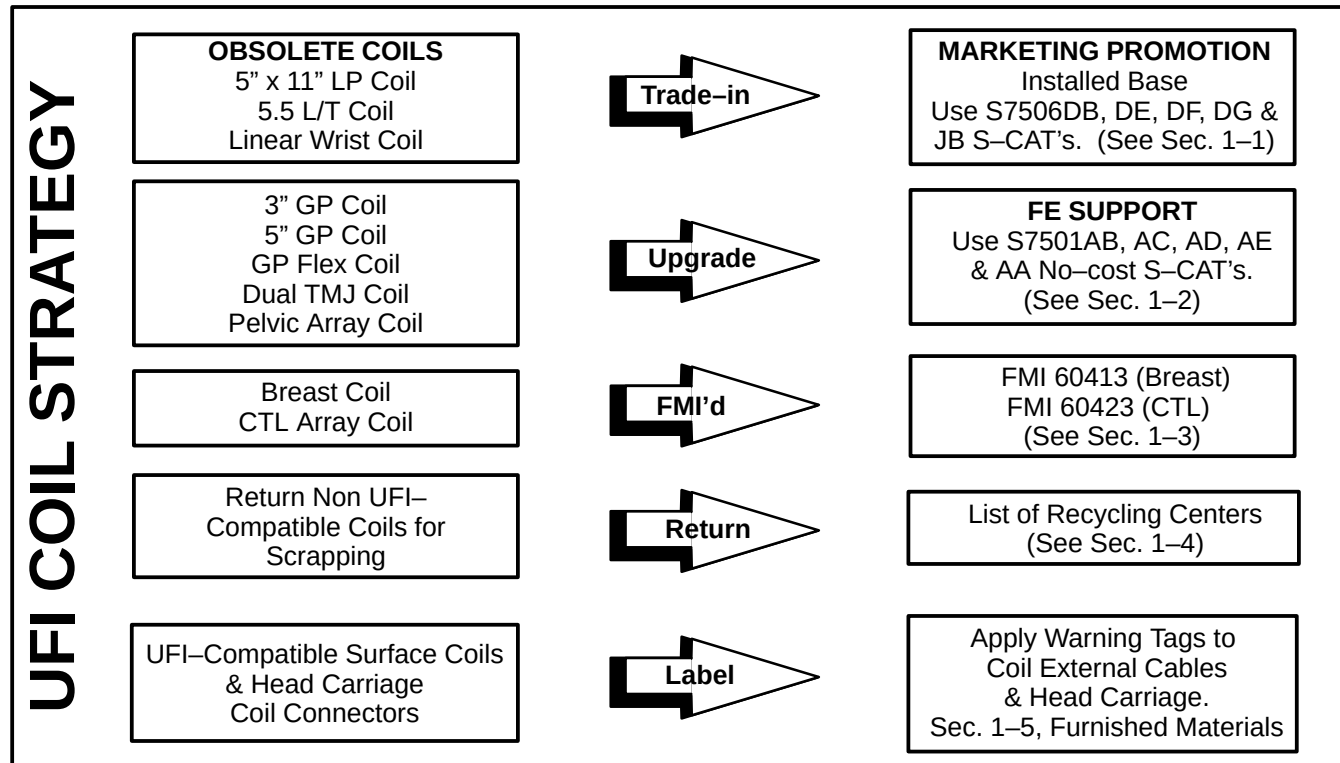
Note

These procedures **do not** apply to 1.5T Signa "Advantage" 5.5 Software Upgrade only systems. For these systems, selecting the "Advantage" Gradient Configuration file option during Load from Cold, will set B1 RMS to the appropriate value.

For NON-GE Approved surface coils: Refer the customer to the original coil manufacturer for information on coil compatibility with Signa Horizon.

Note

GE makes no claims as to compatibility of coils not approved by GE for use on Signa systems. We have no way of knowing if any of these non-GE coils are compatible with 5.5 software and/or hardware. The customer uses these coils at their own risk, and GE accepts no responsibility for their use.

Service Readiness:**1-1 MARKETING PROMOTION**

The Following Coils Are Not UFI-compatible. These coils should **NOT** have any UFI-compatible warning tags attached: (See Table.)

The 5" x 11" LP, 5.5" L/T, and Linear Wrist Coils will be replaced through an Installed Base Marketing Promotion to the customer. The customer will receive a letter that will allow him to trade-in the LP and/or L/T coils for a GP Flex coil or Quad T/L coil. The customer will be able to trade-in the Linear Wrist coil for a Quad Wrist coil (E-CAT). This promotion must occur through the GE Sales force. Trade-in allowances and prices are in effect for orders placed through December 31, 1997 and installed by March 31, 1998. This promotion expires on December 31, 1997.

Each of these Trade-in Allowances requires the physical return of a customer owned surface coil(s). Depending on whether one or two coils are being traded-in, the GE field engineer will remove the existing non-compatible coil(s) at the time of the new coil installation.

UFI Non-Compatible Coils:

The non-compatible coils should be returned to your local recycling center for scrapping (See Section 1-4).

COIL (1.5T)	OLD M-CAT# (DELISTED)	OLD GE #
1.5T 5" X 11" LP (Obsolete)	M1085AF	46-258882G1 46-306773G1
1.5T 5.5" Planar L/T (Obsolete)	M1085AJ	46-265503G1
1.5T Linear Wrist Coil (Obsolete)	M1085BA	46-317224P1 46-328444P1

The FRU strategy for the LP and L/T coils left in the field after the promotion is to replace these coils with a new GP Flex coil. The FRU strategy for the Linear Wrist coils left in the field after the promotion is to replace this coil with a similar Linear Wrist coil. The FRU strategy for the Quad Wrist coils in the field is to replace this coil with a similar Quad Wrist coil.

UFI Compatible Coils:

UFI COMPATIBLE TRADE-IN OPTIONS
With one coil, LP or L/T: Upgrade to GP Flex; Order S7506DB Upgrade to Quad T/L; Order S7506DF
With both LP and L/T: Upgrade to GP Flex; Order S7506DE Upgrade to Quad T/L; Order S7506DG
Upgrade to Quad Wrist; Order S7506JB

1-2 FIELD UPGRADES

The Following Coils Are Not UFI-compatible. These coils should **NOT** have any UFI-compatible warning tags attached: (See Table on next page.)

The 3" GP, 5" GP, GP Flex, Dual TMJ Array (with integrated 3" GP coils), and Pelvic Phased Array coils can be replaced with new UFI-compatible coils (S-CAT's). All of these affected coils for each upgraded 1.5T Signa Horizon system needs to be identified before or at the beginning of the system installation. Coils shared between other systems at the site are covered by this process, but there is a limit of one (1) of any Pelvic or Dual array coil and two (2) of any 3" GP, 5" GP, or GP Flex used at the site unless further approval is obtained. If the customer wishes to upgrade additional site coils to UFI-compatible coils, the new UFI-compatible coil M-CAT must be ordered through the site sales representative.

For US Customers only: When these coils are identified, contact Bob Healy at GE in Milwaukee and let him know which coils they are (S-CAT's) and how many of each. Bob will add replacement coils onto the customer's FDO for appropriate delivery to the site. **Bob can be contacted by FAX at 414-827-3348, E-Mail HEALY-R.** Bob can also be contacted by phone with any questions at (414) 827-3284.

For Non US Customers: When these coils are identified, contact your site sales representative to verify that these coils (S-CAT's) have been added to the FDO.

UFI Non-Compatible Coils:

The non-compatible coils should be returned to your local recycling center for scrapping (See Section 1-4).

COIL (1.5T)	OLD M-CAT# (DELISTED)	OLD GE #
1.5T 3" Round GP	M1085AG	46-265311G1 46-307144G1
1.5T 5" Round GP	M1085AH	46-258773G1 46-307080G1
1.5T GP Flex Coil	M1085GP	46-317492G1
1.5T Dual TMJ	M1087AD	46-317395G1
1.5T Pelvic Array Coil	M1087AC (FW50 '95)	46-317000G1 46-320300G1 46-328072G1

Note

On M1087AC with M1087AT, the 1.5T ATD-III, remove the interface board, 46-328131P1 or 46-328131P2, from the Pelvic Array for use in your replacement M1087PA. The M1087PA does not contain an interface board.

UFI Compatible Coils, S-CATS:

COIL (1.5T)	NEW M-CAT#	ORDER S-CAT# (UFI GE #)
1.5T 3" Round GP	M1085GA	S7501AB (2127315)
1.5T 5" Round GP	M1085HA	S7501AC (2127316)
1.5T GP Flex Coil	M1085GR	S7501AD (2128554)
1.5T Dual TMJ	M1087DA (Adapter Only)	S7501AE (2127315 (2)) (Adapter + 3" GP coils (2))
1.5T Pelvic Array Coil	M1087PA	S7501AA (2144100)

1-3 FMI's

The Following Coils Are Not UFI-compatible. These coils should **NOT** have any UFI-compatible warning tags attached: (See Table.)

The Breast Coil is being upgraded by a Mandatory FMI. The CTL coil is being upgraded by a Mandatory Safety FMI. Please verify that your site has had both FMI's completed. If not, please contact your local FMI coordinator to request the FMI kits.

UFI Non-Compatible Coils:

The non-compatible coils should be returned in the FMI shipping carton to GPO for upgrading and return to FMI pool **as soon as possible**.

COIL (1.5T)	OLD M-CAT# (DELISTED)	OLD GE #
1.5T Breast Coil	–	46-320785P1
1.5T CTL Coil	M1087AB	46-301800G1 46-320971G1

UFI Compatible Coils:

COIL (1.5T)	NEW M-CAT#	NEW UFI GE#	COIL STATUS
1.5T Breast Coil	M1085BR	46-320785P2	Coil is upgraded by FMI#60413. (In Process) (Mandatory FMI)
1.5T CTL Coil	M1087CT	2132617 2104799	Coil is upgraded by FMI#60423. (Available FW44 '95) (Safety FMI)

1-4 RETURN NON UFI-COMPATIBLE COILS FOR SCRAPPING:**For GEMS-AM:**

Return surface coil(s) to: GE Recycle Center
2200 E. College Ave.
Cudahy, WI. 53110
Att: Paul Neumiller (414) 747-6997

For GEMS-E:

Return surface coil(s) to: Site Recyclage GEMS-E
300 Rue Fourny
78530 Buc – France
Att: Annie Thomas, Phone # 30704742

For GEMS-A:

Return surface coil(s) to: Your country's local GEMS Recycling Center or one of
the above
recycling centers.

1-5 FURNISHED MATERIALS

A kit, 2146243 (FRU-2), is attached to inside front frame of the RF/PEN Cabinet above the RF System Controller Module. This kit includes:

Item	Part Number	FRU	Name	Quantity	Description (Remarks)
1	2145936	X	INSTALLATION	1	INSTALLATION DIRECTIONS
2	2146833	2	VINYL STRIPS	20	WHITE VINYL STRIP, 1" X 5", ADHESIVE BACK
3	2146832	2	WARNING TAG	20	UFI WARNING TAG, ADHESIVE BACK
4	2146831	2	WARNING LABEL	2	UFI WARNING LABEL, ADHESIVE BACK, ENGLISH
5	2146831-2	2	WARNING LABEL	2	UFI WARNING LABEL, ADHESIVE BACK, SPANISH
6	2146831-3	2	WARNING LABEL	2	UFI WARNING LABEL, ADHESIVE BACK, PORTUGUESE
7	2146831-4	2	WARNING LABEL	2	UFI WARNING LABEL, ADHESIVE BACK, ITALIAN
8	2146831-5	2	WARNING LABEL	2	UFI WARNING LABEL, ADHESIVE BACK, GERMAN
9	2146831-6	2	WARNING LABEL	2	UFI WARNING LABEL, ADHESIVE BACK, FRENCH
10	2146831-7	2	WARNING LABEL	2	UFI WARNING LABEL, ADHESIVE BACK, CHINESE
11	2146831-8	2	WARNING LABEL	2	UFI WARNING LABEL, ADHESIVE BACK, JAPANESE