



GE Medical Systems

Technical Publications

Direction 2320276
Revision 0

Signa® LX Intel 8 Port Fast Ethernet Hub Upgrade

Copyright© 2002 by General Electric Co.

Operating Documentation

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation "**damage in shipment**" written on **all** copies of the freight or express bill **before** delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

Immediately complete a "Damage Loss Claim Form", available via MS Exchange Mail, after the damage is found.

MS Exchange Path:

Outlook/Public Folder/All Public Folders/Medical Systems/!Global Initiatives/Information Management/Forms/Common Forms/DAMAGE LOSS CLAIM FORM.

Send the completed form to the email address listed in the form.

For more information about the Transportation Claim Procedure, access the GE Medical Systems Intranet and enter the following URL address (case sensitive):

<ftp://3.87.40.2/globepro/qualsys/Docs/190016MF.PDF>

Rev. 11/15/2000

REVISION HISTORY

<u>REV</u>	<u>DATE</u>	<u>PRIMARY REASON FOR CHANGE</u>
------------	-------------	----------------------------------

0	Jan 30, 2002	Initial version for Release.
---------	--------------	------------------------------

<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>
-------------	------------	-------------	------------	-------------	------------	-------------	------------

Title Page	0						
------------------	---	--	--	--	--	--	--

2-6.....	0						
----------	---	--	--	--	--	--	--

* This revision/letter corresponds to the indicated document's revision control system.

Blank page

This Direction documents the upgrade installation of Intel 10/100 Base T Fast Ethernet Hub Kit 2272740-2 on an Octane Host Computer workstation, which replaces the end of life Ethernet Hub 2270639.

Table 1 provides the list of what is included in the 2272740-2 10/100 Base T Intel 8 port hub kit.

TABLE 1
2272740-2 10/100 BASE T INTEL 8 PORT HUB KIT

Part Number	Quantity	Description
2282293	1	Intel 8 Port 10/100 Fast Ethernet Hub
2310044	1	MR LX Hub Panel With Foam
2320276	1	Installation Instructions

1. Log out from Host Computer.
2. Shut down the Host Computer
3. Bring the PC down by closing the Program Manager.
4. As applicable, open the circuit breaker(s) that power the Operator Workspace to turn power off.



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE OPERATOR WORKSPACE CABINET BEFORE CONTINUING THIS PROCEDURE.

5. Lock and tag the circuit breaker(s). (Refer to CD-ROM *Dir. 2160623 MR Release Signa Lx Service Methods, Renewal Parts and Service Tools, Safety*, Section 6, OSHA LOCKOUT/TAGOUT REQUIREMENTS.)
6. Open computer cabinet door and verify that power to the Operator Workstation is off by twice pressing (on/off) the power button on the Host Computer; the power LED should remain unlighted.
7. Remove the covers from the cabinet (Refer to 2160623 Signa LX Service Methods for guidance's in removing covers). Dismantle the power plug retainers and power strip used to power the Ethernet hub. See Illustration 1. You will be discarding the retainers and bracket but power strip is re-installed.

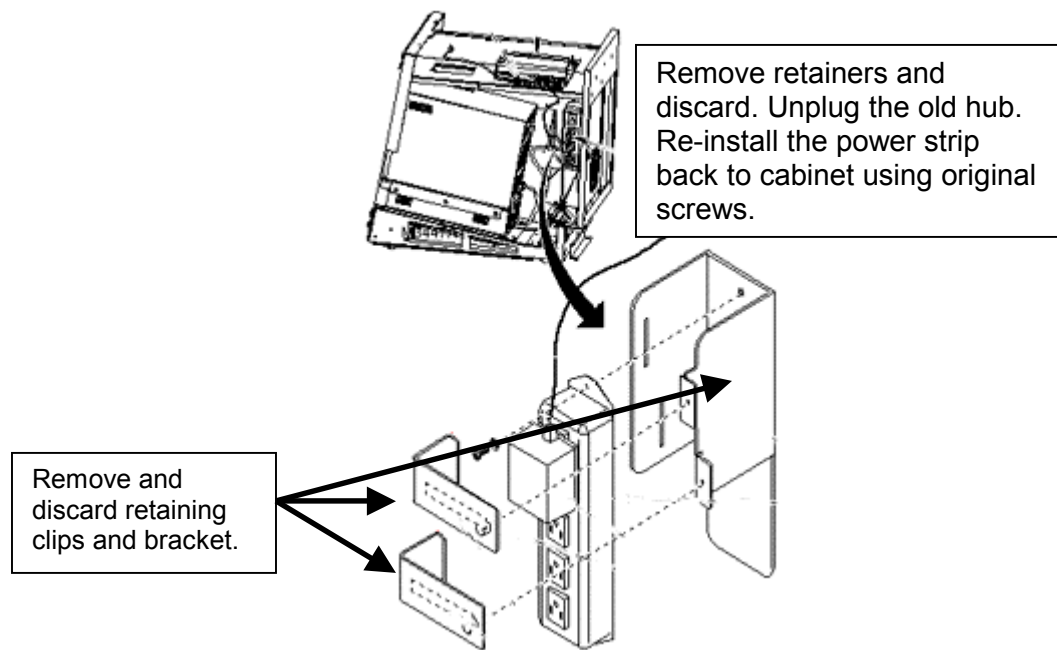


ILLUSTRATION 1
REMOVE UNNEEDED HARDWARE

8. Remove and discard the Hub and hub-restraining bracket holding it in place. See Illustration 2.

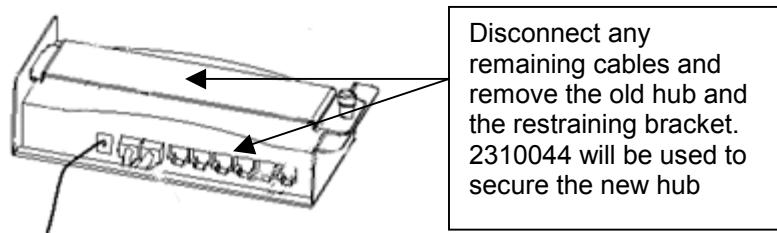


ILLUSTRATION 2
REMOVE HUB

9. Install new hub and secure with new 2310044 panel with foam. Connect the power plug for the hub to the power strip. Reconnect run 811 from the Octane computer to port 2 and the system PC to port 6.
10. Reassemble the WorkSpace cabinet, apply power and test system for proper operation.