

TABLE OF CONTENTS

TABLE OF CONTENTS 1

1- LASER LIGHT CHECKS 2

REVISION HISTORY 5

1- LASER LIGHT CHECKS

For systems with Non-Cx magnet refer to procedure for Patient Alignment Light Checks for Non-Cx Magnet.

Tools Required


- SPT Short Loader Cylinder, 2135652

Note

SPT Cylinder Availability - The SPT Short Loader Cylinder is part of the SPT Short Loader Assembly (2125244). This hardware, though shipped with Signa systems, is the property of GE. It is available for customer purchase.


The three laser lights, mounted on the front face of the magnet bore, are preadjusted at the factory and do not normally require realignment after the initial installation alignment procedures are completed.

1. Place SPT Short Loader Cylinder on patient table.
2. Press Align On button on magnet enclosure.
3. In the event that the laser beams are not aligned as shown in Illustration L2603A and Illustration L2604B, refer to procedure for Patient Alignment Lights, Section 1- Laser Light Alignment Procedures.

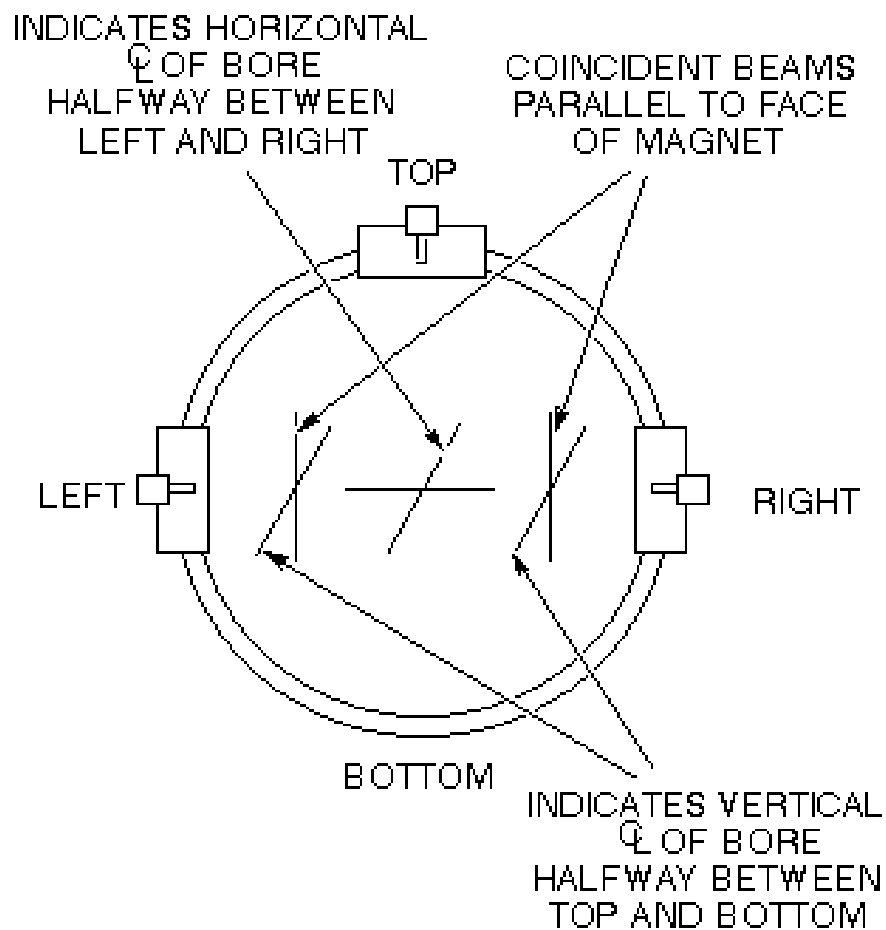

L2603A

HALOGEN LIGHT BEAM PATTERN
ILLUSTRATION L2603A

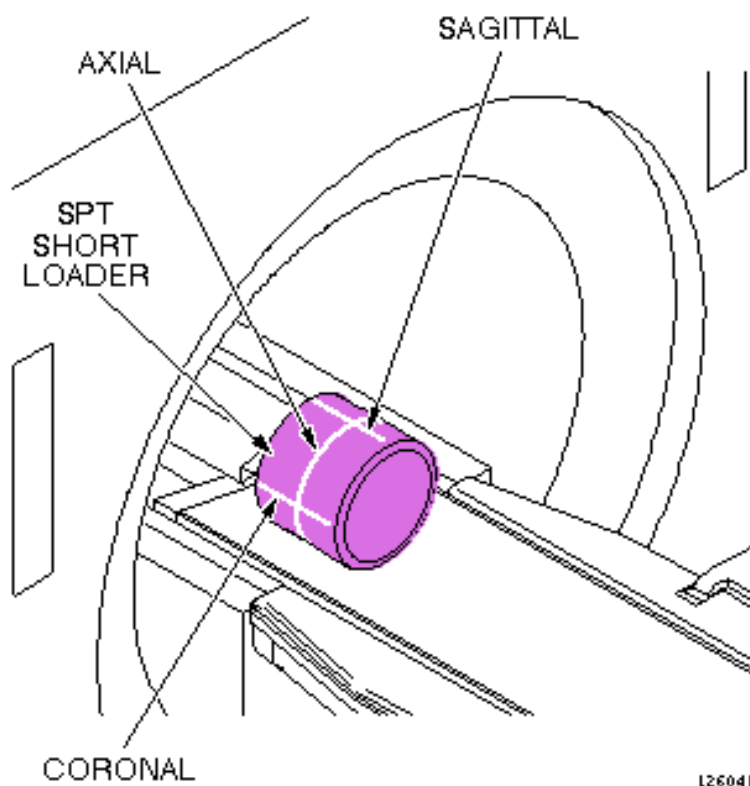

HALOGEN WRAP-AROUND BEAM PATTERN
ILLUSTRATION L2604B

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	June 15, 1998	J. Saperstein	Initial conversion from Toolbook to Word formats.