

TABLE OF CONTENTS

TABLE OF CONTENTS 1

1- ENCLOSURE LOCATIONS/CONNECTIONS 2


REVISION HISTORY 5

1- ENCLOSURE LOCATIONS/CONNECTIONS

For front view of magnet enclosure with component locations, see
Illustration L2692B - Non-Cx Magnets

MG2 COMPONENT DESIGNATORS

- A4** Operator Button Panel & Display (left)
- A5** Operator Button Panel & Display (right)
- A11** Physical Acquisition Controller (PAC)
- A29** Dock Assembly
- A30** Patient Alignment Lights
- A31** Dock Limit Module
- A33** Scan Room Interface Module (SR)
- A36** Bore Temperature Sensor

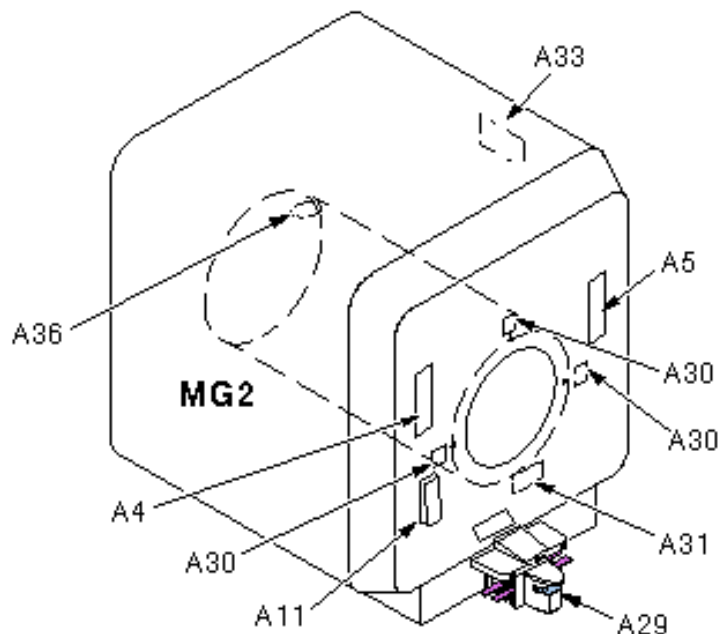


L2692B

FRONT VIEW OF MAGNET ENCLOSURE WITH COMPONENT LOCATIONS (SIGNA HORIZON)
ILLUSTRATION L2692B

Illustration L2692C - Cx Magnet

MG2 COMPONENT DESIGNATORS

- A4** Operator Button Panel & Display (left)
- A5** Operator Button Panel & Display (right)
- A11** Physical Acquisition Controller (PAC)
- A29** Dock Assembly
- A30** Patient Alignment Lights
- A31** Dock Limit Module
- A33** Scan Room Interface Module (SR)
- A36** Bore Temperature Sensor

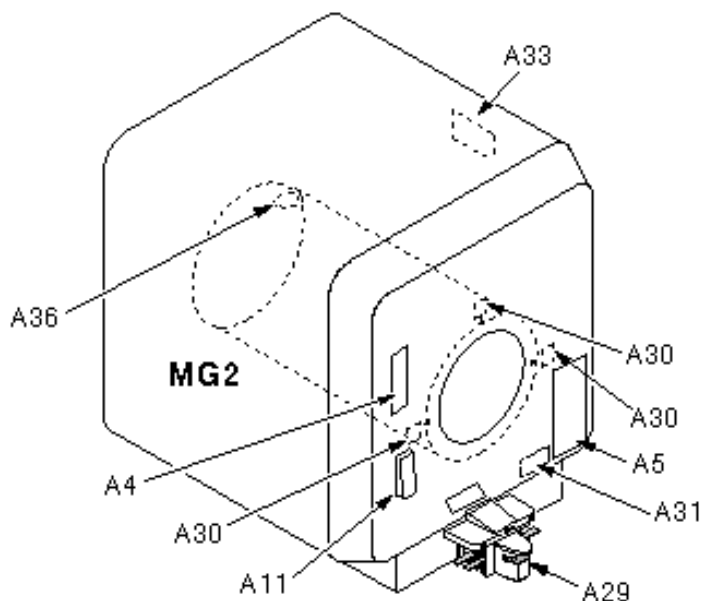

L2692C

FRONT VIEW OF CX MAGNET ENCLOSURE WITH COMPONENT LOCATIONS (SIGNA HORIZON)
ILLUSTRATION L2692C

For rear view of magnet enclosure with component locations, see

Illustration L2693B - Non Cx Magnets

Illustration L2693C - Cx Magnet


For magnet enclosure interconnects, refer to *Direction 2117724, Signa 5.x Block Diagrams and Supplemental Schematics*, Patient Handling tab.

MG3 COMPONENT DESIGNATORS

- A1** Magnet Interface #1
A1 Limit Switch Board
- A2** Patient Communication Box
A1 Patient Communications Board
- A6** Longitudinal Encoder
- A7** Longitudinal Motor Drive
- A12** Longitudinal End of Travel
(Only on some Upgrades)
- A17** Head Coil Cable Takeup Interface

MG2 DESIGNATORS

- A14** Patient Comfort Module
- A15** Bore Light
- A16** Head Coil Carriage Assembly
 - A1 Head Coil (NOT SHOWN)
 - A3 Head Preamplifier
 - A4 Quick Disconnect
 - A5 Head Coil T/R Switch
 - A6 Coil Adapter
- A18** Bias Bulkhead

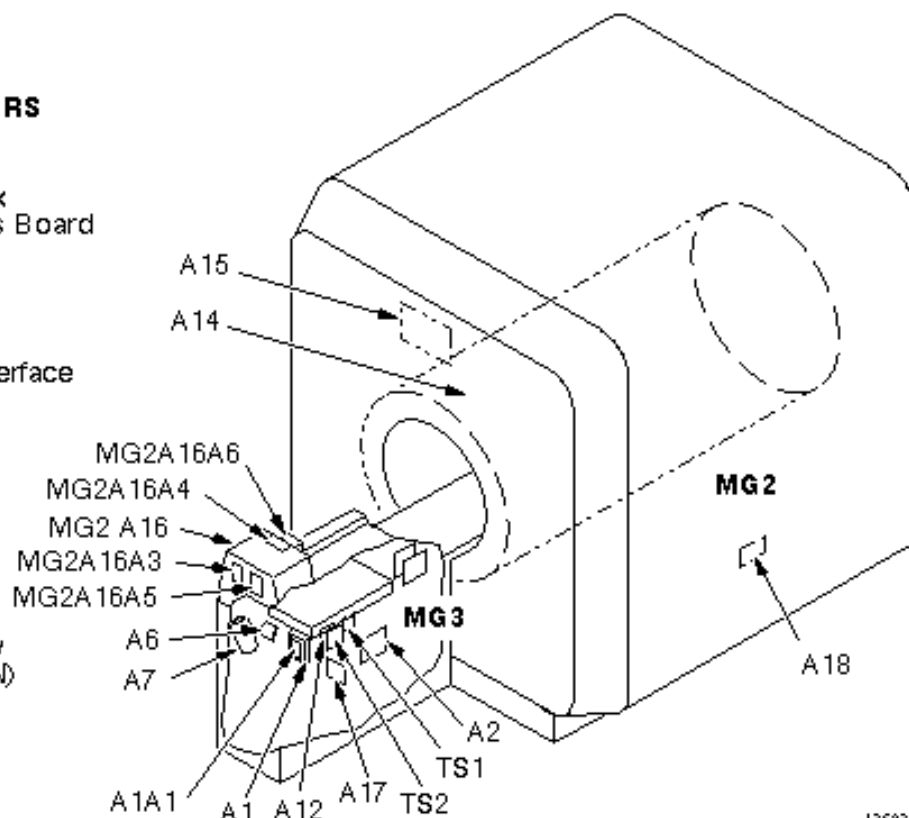

L2693B


REAR VIEW OF MAGNET ENCLOSURE WITH COMPONENT LOCATIONS (SIGNA HORIZON)
ILLUSTRATION L2693B

MG3 COMPONENT DESIGNATORS

- A1** Magnet Interface #1
A1 Limit Switch Board
- A2** Patient Communication Box
A1 Patient Communications Board
- A6** Longitudinal Encoder
- A7** Longitudinal Motor Drive
- A12** Longitudinal End of Travel
(Only on some Upgrades)
- A17** Head Coil Cable Takeup Interface

MG2 DESIGNATORS

- A14** Patient Comfort Module
- A15** Bore Light
- A16** Head Coil Carriage Assembly
 - A1 Head Coil (NOT SHOWN)
 - A3 Head Preamplifier
 - A4 Quick Disconnect
 - A5 Head Coil T/R Switch
 - A6 Coil Adapter
- A18** Bias Bulkhead

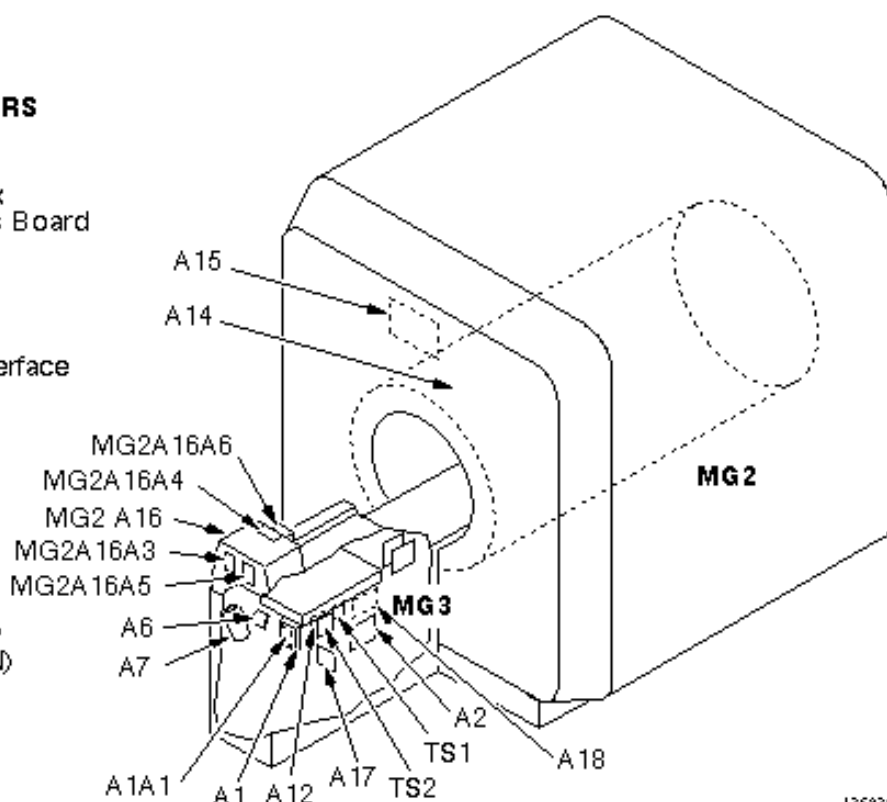

L2693C

REAR VIEW OF CX MAGNET ENCLOSURE WITH COMPONENT LOCATIONS (SIGNA HORIZON)

ILLUSTRATION L2693C

For connections to/from the magnet enclosure, you can also refer to the system interconnect map (provided with your system installation manual).

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	June 15, 1998	J. Saperstein	Initial conversion from Toolbook to Word formats.