

TABLE OF CONTENTS

TABLE OF CONTENTS 1

1- PAC2 BOARD SETUP 2

2- PAC2 BOARD ADJUSTMENTS REQUIRED 4

3- PAC2 BOARD FUNCTIONAL CHECKS REQUIRED 4

REVISION HISTORY 6

Description - Procedure for replacing the Physiological Acquisition Controller (PAC2) board on Signa Horizon systems.

1- PAC2 BOARD SETUP



Equipment damage possibility. Before replacing the PAC2 assembly, disconnect power to avoid damaging the PAC2 board.

1. Lock out and tag out circuit breaker MR1 A14 CB6 at rear of RF/Pen Cabinet on magnet enclosure power supply (MEPS) module.

For **RF/Pen II** and **RF/PDU Cabinet**: Lock out and tag out circuit breaker MR1 A20 CB4 at rear of RF/Pen II and RF/PDU Cabinet (SSM).

(Refer to CD-ROM *Dir. 2187583-3 [or -2], MR Release Signa 5x/8x Service Methods, Renewal Parts and Service Tools*, Procedure for Safety, Section 6, OSHA LOCKOUT/TAGOUT REQUIREMENTS.)

2. Verify that voltage between MR1 A14 TP4 (+) and TP13 (–) is 0 Vdc.

For **RF/Pen II** and **RF/PDU Cabinet**, verify that voltage between MR1 A20 TP4 (+) (SRI) and TP13 (–) (SRI RTN) is 0 Vdc.

3. Remove PAC2 assembly from enclosure.

4. For board information, refer to:

- Illustration L2928B - Board layout
- Table 1 - Switches
- Table 2 - LEDs
- Table 3 - Test points

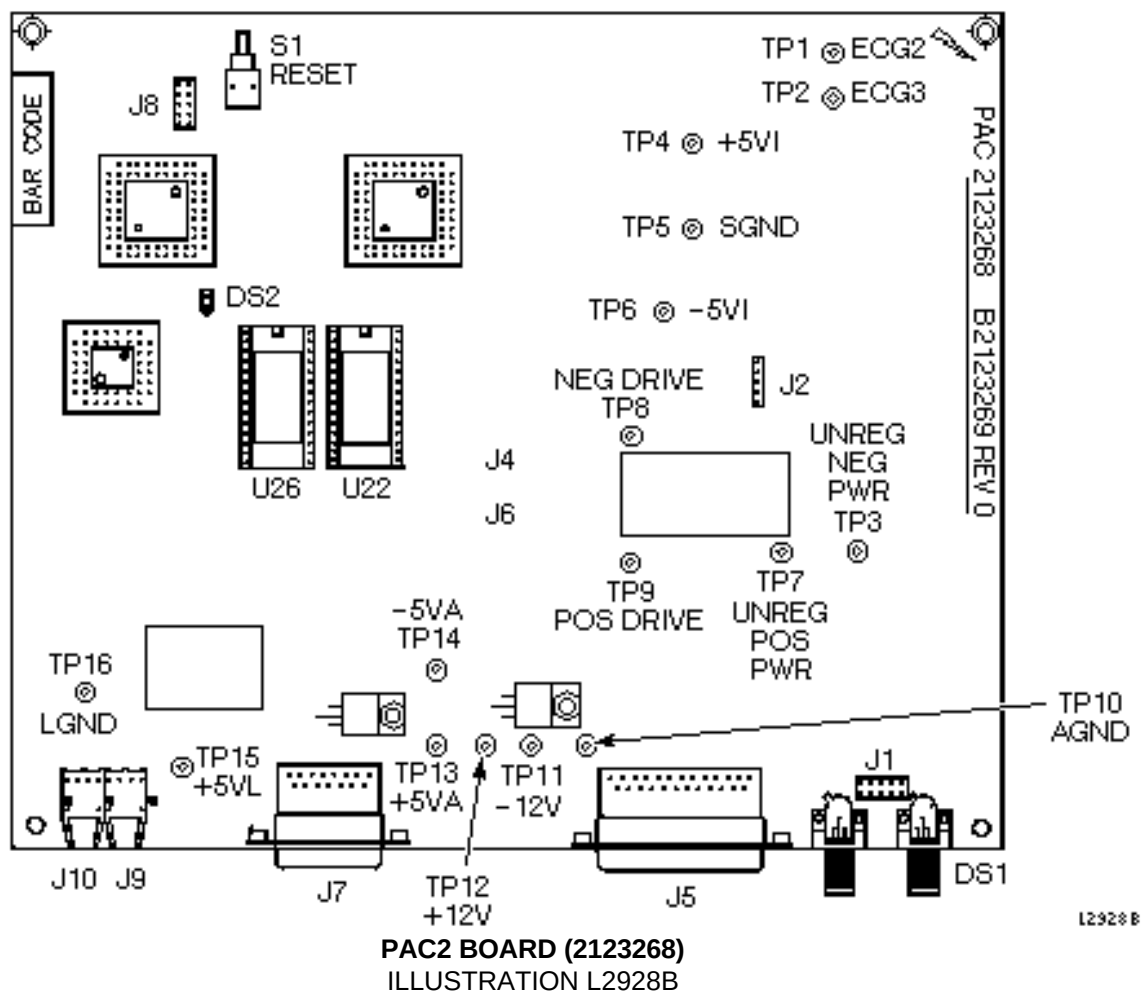


TABLE 1
SWITCH INFORMATION

Switch	Position	Description	Notes
SW1	Momentary	Reset	Reset power-up switch

TABLE 2
LED INFORMATION

LED	Color	Description	Notes
----		Bar graph	Respirational signal
DS2	Green	Operation/ Diagnostics	
----		Bar graph	Peripheral pulse signal

TABLE 3
TEST POINT INFORMATION

Test Point	Name	Notes
TP1	ECG2	ECG 2 Lead input to ISO-A/D
TP2	ECG3	ECG 3 Lead input to ISO-A/D
TP3	UNREG NEG PWR	Unregulated negative ISO-power (+5 Vdc $\pm$ 10%)
TP4	+5VI	Isolated +5V supply
TP5	SGND	Isolated power ground
TP6	+5VI	Isolated -5V supply
TP7	UNREG POSPWR	Unregulated positive ISO-power (+10 Vdc $\pm$ 10%)
TP8	NEG DRIVE	Air core transformer negative drive
TP9	POS DRIVE	Air core transformer positive drive
TP10	AGND	Analog ground
TP11	-12V	-12 Volt analog supply
TP12	+12V	+12 Volt analog supply
TP13	+5A	+5 Volt supply
TP14	-5A	-5 Volt supply
TP15	+5L	+5V VCC Power bus
TP16	LGND	Logic ground

5. Install replacement PAC2 assembly.



Equipment damage possibility. Do not connect power to the PAC2 assembly until the calibration procedure has been completed. Connecting power before calibration may damage the PAC board.

2- PAC2 BOARD ADJUSTMENTS REQUIRED

1. Refer to procedure for PAC Power Supply Calibration. Perform all procedures.
2. Remove tag-out tag from circuit breaker MR1 A14 CB6 at rear of RF/Pen Cabinet on MEPS module.
For **RF/Pen II** and **RF/PDU Cabinet**: Remove tag-out tag from circuit breaker MR1 A20 CB4 at rear of RF/Pen II and RF/PDU Cabinet (SSM).
3. Restore power to PAC2 by switching on circuit breaker MR1 A14 CB6, or MR1 A20 CB4 at rear of **RF/Pen II** and **RF/PDU Cabinet** (SSM).

3- PAC2 BOARD FUNCTIONAL CHECKS REQUIRED

1. Refer to procedure for Respiratory Comp & Bellows Check and procedure for ECG Leads Check, respectively. Perform all procedures.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	June 16, 1998	J. Saperstein	Initial conversion from Toolbook to Word.