

TABLE OF CONTENTS

TABLE OF CONTENTS.....	1
1- +12VDC BORE LIGHT POWER SUPPLY	2
2 - PATIENT INTERCOM MICROPHONE BIAS ADJUSTMENT	2
REVISION HISTORY.....	4

Description - Procedures for adjusting the Bore Light Power Supply and microphone bias voltage.


1- +12VDC BORE LIGHT POWER SUPPLY

RF/Pen Cabinet: The Bore Light Power Supply is located on the High Current Power Supply Board (MR1 A14 A2) in the RF/Pen Cabinet. To adjust the Bore Light Power Supply voltage, refer to procedure for Magnet Enclosure P.S. Module Calibration.

RF/Pen II and RF/PDU Cabinets: The Bore Light Power Supply is located on the Dock/Light Board (MR1 A20 A5 A2) in the RF/Pen II and RF/PDU Cabinets. There are no bore light adjustments that can be done.

2 - PATIENT INTERCOM MICROPHONE BIAS ADJUSTMENT

The bias voltage for the patient microphone must be adjusted to obtain proper microphone response. This voltage is adjusted on the Remote Intercom Board (MG3 A2 A1) inside the Patient Communication Box (MG3 A2) on the rear pedestal. See Illustration 2-1.

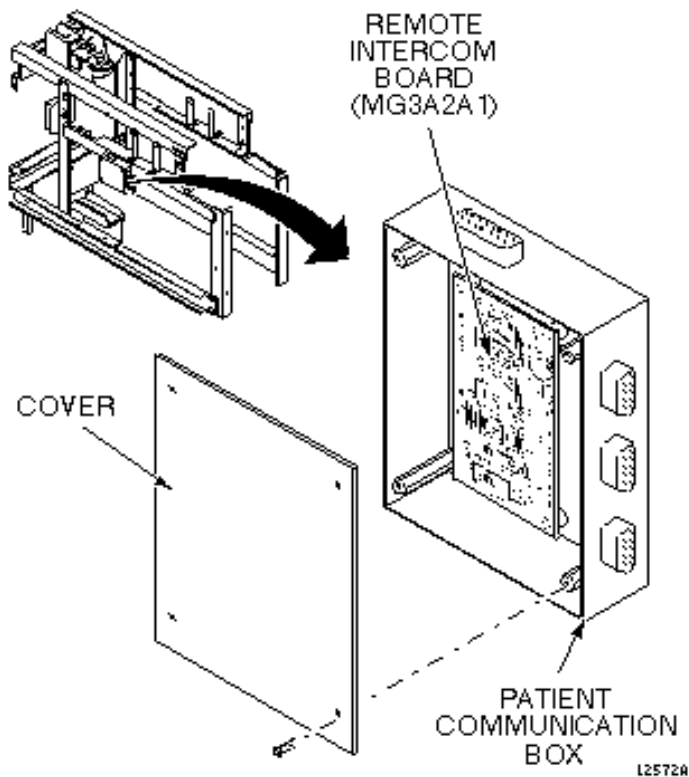


PATIENT COMMUNICATION BOX
ILLUSTRATION 2-1

Use the following procedure to adjust bias voltage:

1. Ensure that bias cable MG3 A2 J6 is connected to the patient communication box.
2. Remove Patient Communication Box (MG3 A2) cover to expose Remote Intercom Board. See Illustration 2-1.

3. Connect voltmeter leads across capacitor C42 on the Remote Intercom Board (MG3 A2 A1), as shown in Illustration 2-2.

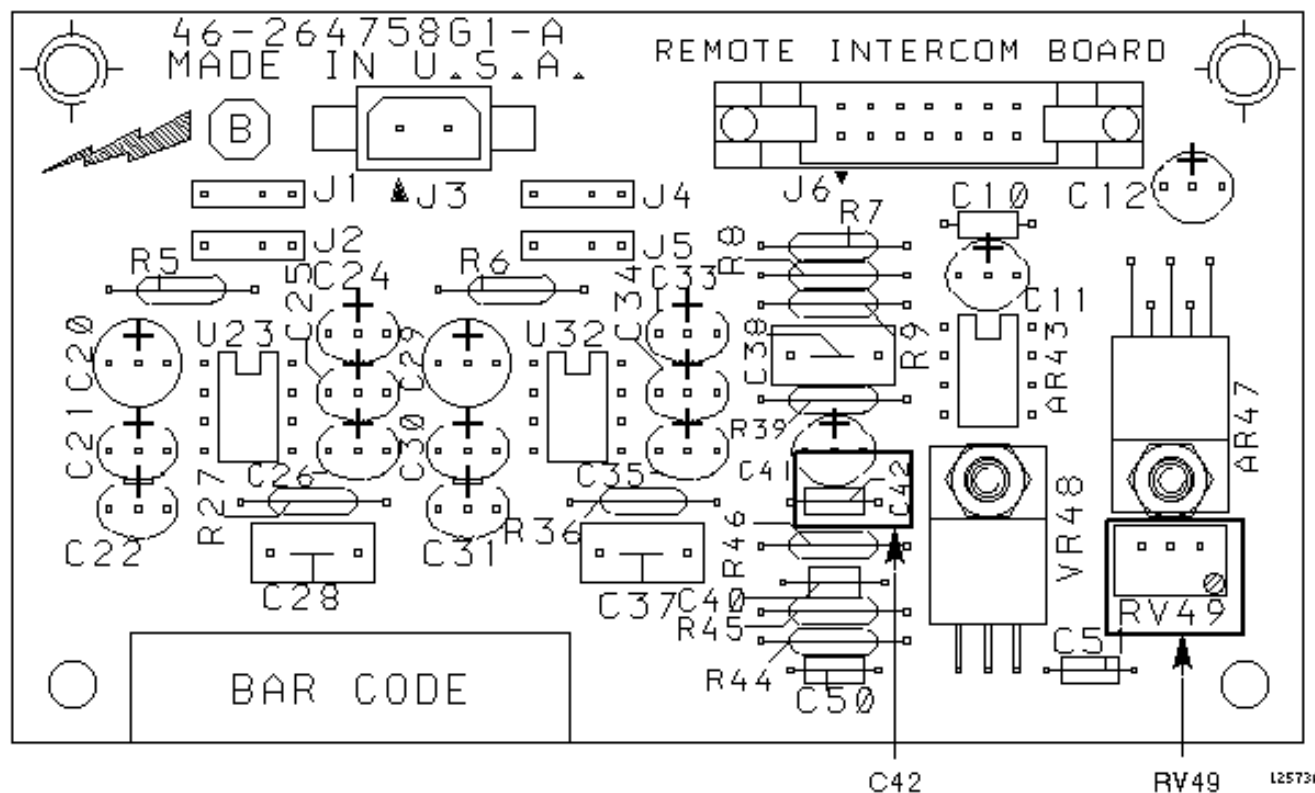

REMOTE INTERCOM BOARD ILLUSTRATION 2-2

4. Adjust RV49 (Illustration 2-2) to obtain voltage of $+6.0 \pm 0.2\text{Vdc}$.
5. Replace Patient Communication Box cover.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	June 16, 1998	J. Saperstein	Initial conversion from Toolbook to Word.
1	January 6, 1999	K. Keshena	Updates to styles.