

TABLE OF CONTENTS

TABLE OF CONTENTS 1

1- SCAN-CONTROL INTERCOM MODULE (SCIM) AND KEYBOARD INSTALLATION..... 2

 1-1 Installing The SCIM and Keyboard..... 2

A- APPENDIX - SCIM BUTTON IDENTIFICATION 5

REVISION HISTORY 6

1- SCAN-CONTROL INTERCOM MODULE (SCIM) AND KEYBOARD INSTALLATION

The Scan-Control Intercom Module (SCIM) is the current keyboard on all Signa LX and OpenSpeed systems. This keyboard is designed to have the same form and function of the keyboard on the CT systems adding uniformity in operation between the MR and CT Operator control functions.



SCAN-CONTROL INTERCOM MODULE (SCIM) AND KEYBOARD
ILLUSTRATION 1-1


1-1 Installing The SCIM and Keyboard

1. Remove the SCIM, Data Cable, and Keyboard From the box.



SCIM KIT CONTENTS
ILLUSTRATION 1-2


2. Flip the SCIM unit up side down exposing the cable connections.

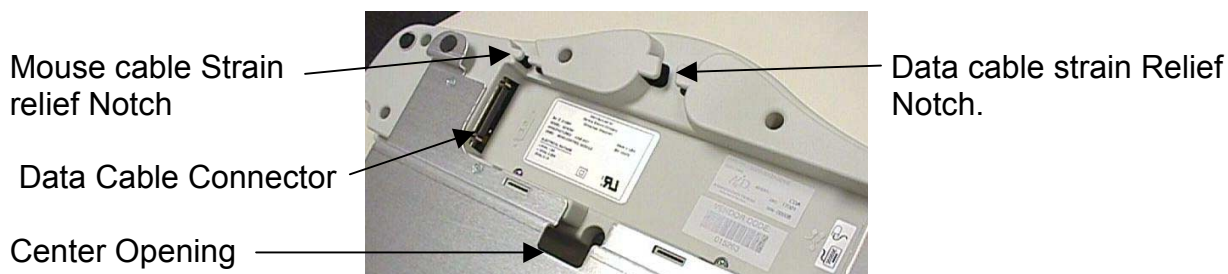

DATA CABLE CONNECTOR
ILLUSTRATION 1-4


3. Attach the following cables:

- a. Attach the data cable to the new SCIM Scan-Control Intercom Module by plugging in the side of the cable that is marked " Keyboard". Route the cable through the strain relief notches in at the rear of the SCIM. (Refer to Illustration 1-4)
- b. Attach the mouse to the jack labeled with a Mouse icon. Route the cable through the strain relief notches in at the rear of the SCIM. (Refer to Illustration 1-5)

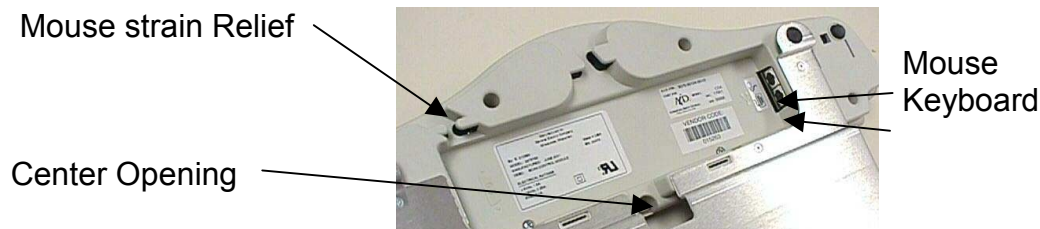

MOUSE AND KEYBOARD CONNECTIONS
ILLUSTRATION 1-5

- c. Flip the SCIM over so it is right side up and guide the keyboard cable through the center opening. (Refer to Illustration 1-5 for location of center opening)
 - d. Once again, flip the SCIM unit up side down and attach the keyboard cable to the jack labeled with a Keyboard icon (Store excess cable in the base of the SCIM by tucking it in the depressed area). (Refer to Illustration 1-5)
4. Flip the SCIM over so it is right side up and place it in front of the monitor.

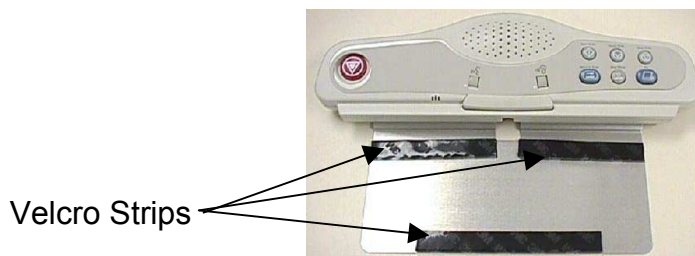

SCIM
ILLUSTRATION 1-6

5. Pull the new keyboard forward, Flipping it on top of the SCIM up side down.
6. Remove the plastic strips from the velcro, exposing the adhesive surface. (Refer to illustration 1-6)
7. Guide the cable back into the center notched opening until all the cable is under the SCIM base.
8. Place the keyboard tight against the SCIM, Center it as best as possible, and drop into place. The Velcro pads will now be sticking to the bottom of the keyboard.

Note


The adhesive on the velcro takes 24 hours to completely adhere to the bottom of the keyboard. Allow the adhesive to dry completely before removing keyboard from SCIM base.

**SCIM AND KEYBOARD ASSEMBLED****ILLUSTRATION 1-7**

9. Position the SCIM Assembly in front of the monitor and guide the data cable through the opening provided in the Workspace Desktop.
10. Connect the data cable to J6 of the Workspace Interface Module (WIM)

This completes the installation of the SCIM Assembly.

A- APPENDIX - SCIM BUTTON IDENTIFICATION

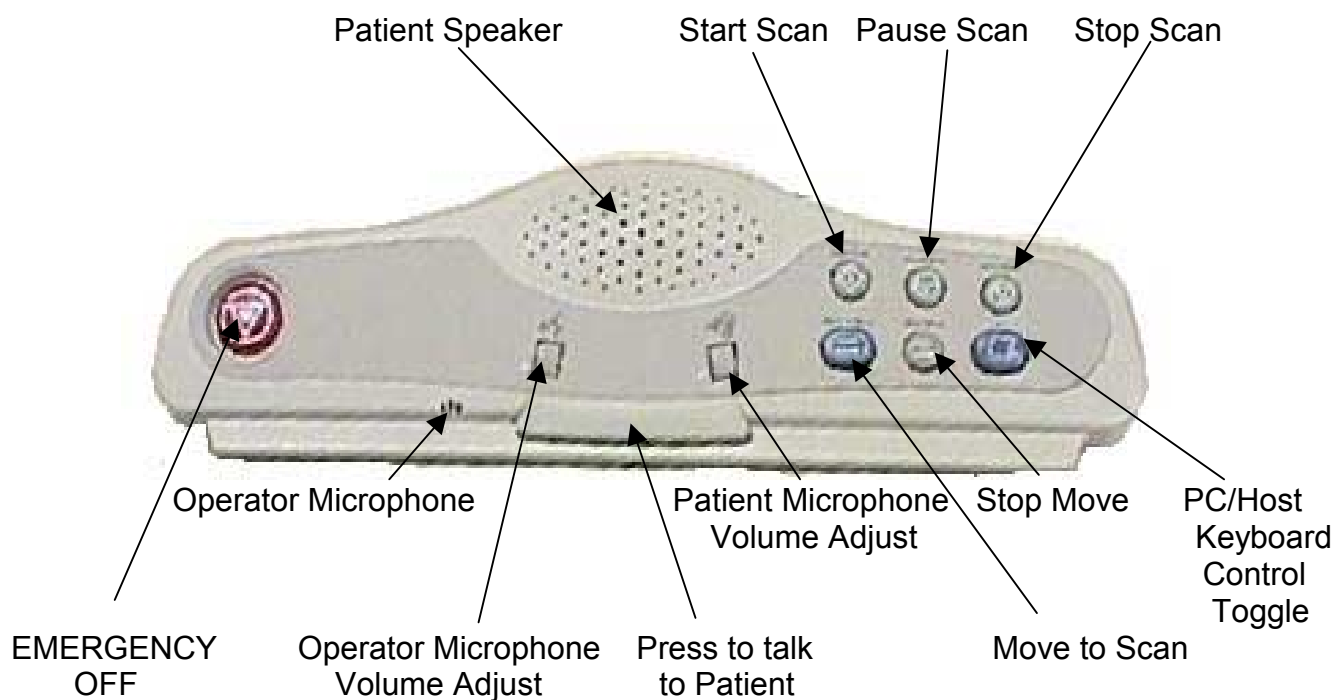

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
A	Oct 3, 2001	D. Hofstetter	Initial version.
1	Dec 19, 2001	D. Hofstetter	Changes per request of LSI.