

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1- INTRODUCTION	3
1-1 Overview.....	3
1-2 Camera Vendor Participation	3
2- NEC 1880SX MONITOR INSTALLATION	3
2-1 Installation of the NEC1880SX in a Fixed Site Environment.....	3
2-2 Installation or Replacement of an NEC 1880SX in a Mobile Environment.....	5
3- VERIFY AMBIENT LIGHTING CONDITIONS.....	6
4- SETTING LCD MONITOR FREQUENCY RATE AND RESOLUTION.....	6
4-1- Setting the Line Rate and Resolution.....	6
5- SIGNA HOST GAMMA SETUP- NEC 1880SX	7
5-1 Gamma Configuration for all SGI Computers	9
5-1-1 Setting Gamma with the Patch CDROM – All Software versions 9.0 and below	9
5-1-2 Setting Gamma - 9.1, HF02 and HFO3.....	10
6- LCD AUTO ADJUSTMENT AND OSM™ ON SCREEN MANAGER OPERATIONS.....	12
6-1 Overview.....	12
6-2 Using the OnScreen Manager.....	13
7- FINAL ADJUSTMENT OF CONTRAST AND BRIGHTNESS	14
7-1 Displaying the SMPTE Pattern	14
7-2 Adjusting for proper Window and Level Setting	16
7-3 Fine Tuning the Contrast and Brightness using the SMPTE Pattern	16
8 - CAMERA CALIBRATION.....	17
8-1 DASM Interpolation Setup.....	17
8-2 Camera Imaging Look-up Table.....	18
8-3 Camera Maximum Optical Density.....	18
8-4 Camera Contrast.....	18
8-5 Anatomical Filming	19
9 - TROUBLESHOOTING GUIDE.....	19
9-1 Explanation of Default Gamma Settings	19
9-2 Complaint: Poor SNR, Dynamic Range Lack of Grayscale.	19
9-3 Complaint: Ghosting or Blurring images or Sheering alphanumeric LCD 's ONLY.	19
9-4 Camera Film has changed after Monitor replacement.....	20
9-5 No Picture.....	21
9-6 Image Persistence.....	21
9-7 Image is unstable, unfocused, or swimming is apparent	21
9-8 LED on monitor is not lit (no green or amber color can be seen).....	21
9-9 Display image is not sized properly	21
9-10 Selected resolution is not displayed properly	22
9-11 Diagnostic Image Quality has been traced to problem with LCD setup.	22
9-12 Display image has a green cast to it	22
9-13 Display image ghosts or gray lines in areas of extreme black & White transitions.	22
9-14 Display image Blurry or Blotchy images Described as poor grayscale.	22
9-15 Customer complaints of poor grayscale, Window or level, lack of definition.	22
9-16 LCD FAQ's:	23
9-17 What is Gamma?	23
9-18 How do I set the correct gamma?	24
9-19 Can I use the gamma command to set gamma?.....	24
9-20 How do I determine if the correct gamma Table or setting is installed?	24
9-21 When I type "gamma" in a c-shell I see a number. Is that the gamma setting?.....	24
9-22 When I type "gamma 1.05" in a c-shell, the LCD changes in appearance. Why?	24
9-23 A gamma value of 1.05 is what we always used in the past can we still use it?	24
9-24 My customer does not like the gamma setting that the Gamma Table provides.	25
9-25 Single value optimization technique.	25
A- APPENDIX A.....	26
A-1 Focus, Screen position, Clock Rate and Phase Adjustments.....	26
A-2 Adjusting Focus Using LX Calibration Tools From Toolchest (Alternate).....	26

A-3 Setting up LCD Display for Contrast/Brightness. Using LX Calibration Tools From Toolchest (Alternate) ...	27
A-4 Monitor Adjustments	28
A-5 Contrast and Brightness Adjustment For NEC1880SX LCD Monitor	28
A-6 Control Descriptions- NEC 1880SX LCD Monitor	28
A-7 Front Panel Controls	29
A-8 On-Screen Controls	29
REVISION HISTORY	32

1- INTRODUCTION

1-1 Overview

This procedure will describe how to adjust the LCD display monitor to match the camera. Once the LCD display monitor is re-calibrated, it is essential to check and possibly re-calibrate the camera.

1-2 Camera Vendor Participation

It is recommended that this procedure is performed with the Camera Vendor field engineer, present, should any camera adjustment be necessary. To optimize customer satisfaction, it is also recommended that you have one of the Customer's filming specialists available for the fine tuning and quality review of the film/LCD display monitor conformance.


2- NEC 1880SX MONITOR INSTALLATION

2-1 Installation of the NEC1880SX in a Fixed Site Environment

1. Shutdown Signa and perform a safety lockout/tagout of the Operator Workspace. See Safety procedure to accomplish this.
2. If necessary, remove the old LCD or CRT from the Desktop.
3. Remove the new LCD monitor from its box and its packing and set the NEC 1880SX LCD color monitor on the workspace desktop.
4. Some items will not be used. These can be disposed of locally.


Discard:

1. DVI cable
2. Standard VGA cable
3. PC CDROM drivers

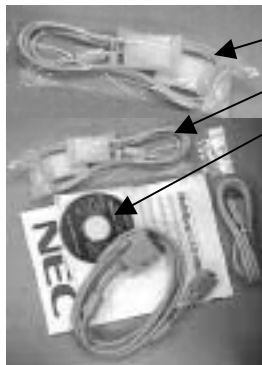

Use:

- Restraint Clips
- Power Cable
- SVGA-RGB Cable

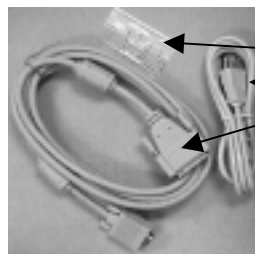

If you are replacing a Host Monitor make sure you also order: GE Part# 2376768
This CDROM disk is labeled "MR Signa® Host Monitor Gamma Correction Software" contains the most up-to date gamma files needed to setup gamma. The software on the CDROM is required, to setup gamma on a replacement monitor. Follow the installation instructions included with the CDROM disk.




GE Gamma Patch
CDROM
(GE Part# 2376768)


DISCARD EXTRA CABLES AND PC CDROM DISK
ILLUSTRATION 2-1

5. Rotate the monitor and connect the Video and power cables.

6. Use the restraint clips provided. Route these cables neatly to the back of the monitor base and the Operators Workspace. See Illustration 2-2 and 2-3.

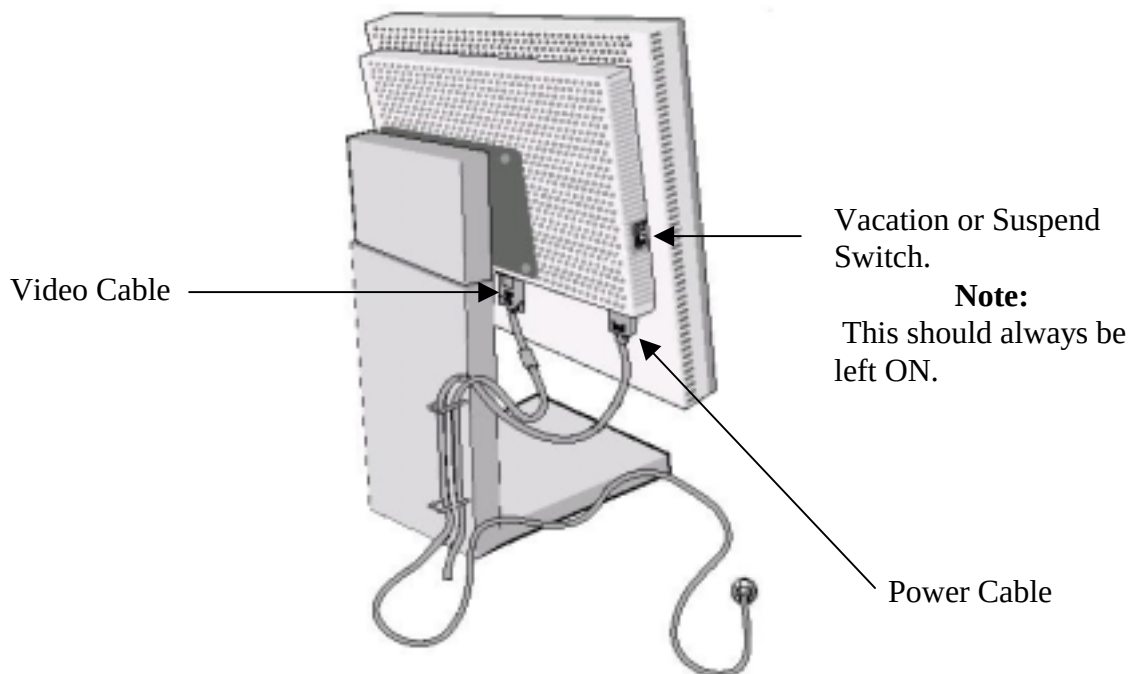

VIDEO AND POWER CABLE CONNECTIONS
ILLUSTRATION 2-2

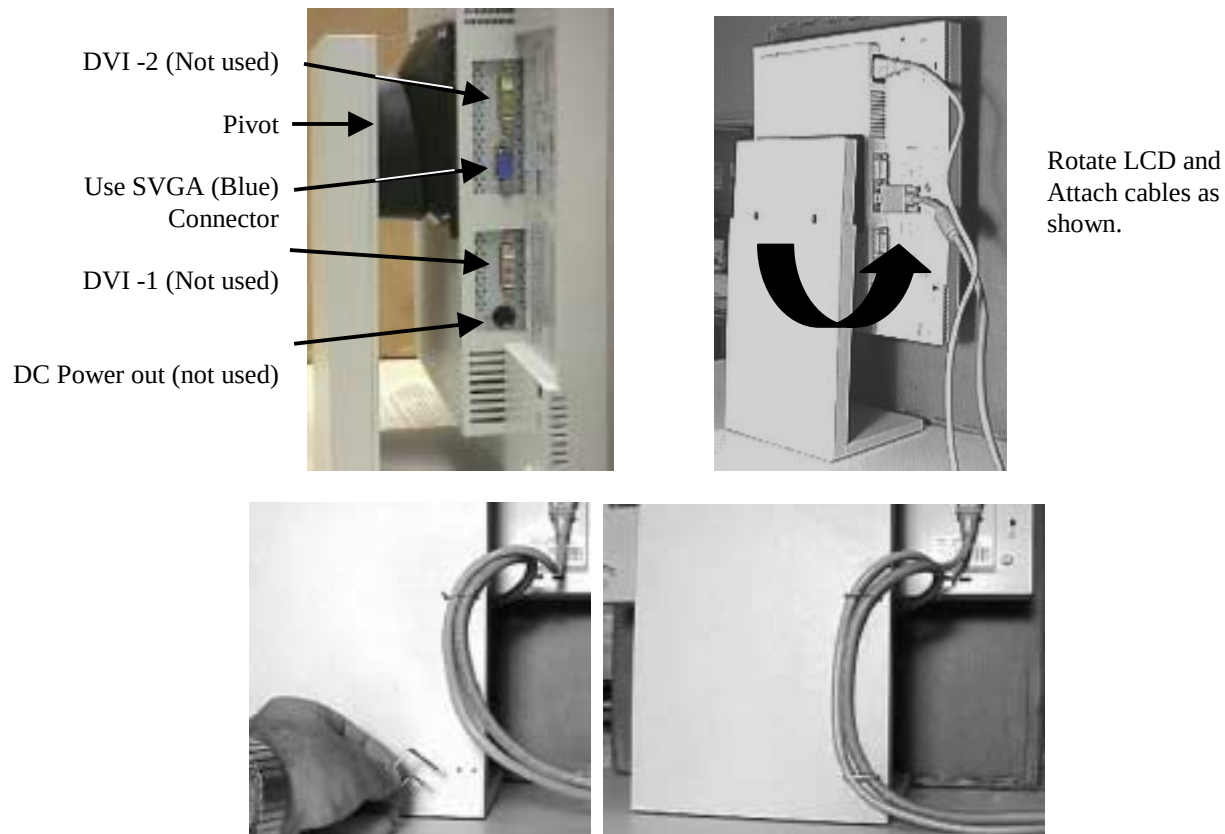

VIDEO AND POWER CONNECTED RESTRAINT CLIPS IN PLACE
ILLUSTRATION 2-3

7. Remove the Safety Lockout of the Operators Workspace and restore power to the system.
8. Verify that main LCD monitor power is ON. If not, press the power button on the front panel and the side panel to turn the main LCD monitor ON. See Illustration 2-3 for side view.
9. The monitor should be positioned no closer than 16 inches and no further away than 28 inches from your eyes. The optimal distance is 24 inches for all monitor types.

Note

Allow the monitor to warm-up for 20 minutes before performing any adjustments.

10. Boot up Signa

11. Password: **adw2.0** <enter>

2-2 Installation or Replacement of an NEC 1880SX in a Mobile Environment

This section pertains to mobile systems only. Replacement of an NEC 1880SX LCD is performed, by removing the four screws that hold the LCD panel to its base. This procedure will only work for the NEC 1880SX.

1. Remove the replacement LCD panel from its shipping box.
2. Place the new LCD Monitor, face down on a flat surface where it will not get scratched or damaged.
3. Remove the 4 screws that hold the LCD panel to its base.
4. Remove the failed NEC 1880SX monitor from its base in the same manner.
5. Using a Phillips screwdriver, mount the new LCD panel to the existing mobile base using the 4 screws.
6. Attach the failed LCD Panel to the new base or dispose of the both items according to normal procedures.

3- VERIFY AMBIENT LIGHTING CONDITIONS

LCD's are extremely susceptible to Viewing Angle. All adjustments and image comparisons should be done from exactly the same viewing angle in all directions. Sit directly in front of the monitor at all times.

In the review area and operator workspace area, verify that the ambient lighting conditions are adjusted to a minimum level. In the operator workspace area, there should be only sufficient light for safely operating the system.

In the review area and operator workspace area, verify that light-boxes are not emitting light, or are properly masked, when not displaying film. This will be a source of excessive glare

In both review area and operator workspace area, verify that there is no source of glare for reviewing films or setting up the images for film. For example, windows should not allow direct light. (blinds should be closed).

Note that both, the operator workspace area and the review area artificial lighting type should be of the same type.

4- SETTING LCD MONITOR FREQUENCY RATE AND RESOLUTION

Failing to do this section will result in Image Ghosting and unnecessary Troubleshooting.

4-1- Setting the Line Rate and Resolution



THIS PROCEDURE NEEDS TO BE DONE EVERY TIME SOFTWARE HAS BEEN RELOADED. LINE RATE AND GAMMA SETTINGS ARE NOT CAPTURED ON THE SAVEINFO MOD AND IS NOT CAPTURED IN THE HEALTH PAGE STATUS OF THE SYSTEM.

MONITOR "LINE RATE" OR "REFRESH RATE", IS REFERRED TO AS "FRAME RATE" AND SWAP RATE", ON SGI SYSTEMS. FOR SIMPLICITY, THE TERM "LINE RATE" WILL BE USED IN THIS DOCUMENT.

Note:

Failure to set the Line Rate to 60Hz for ALL LCD Host monitors will result In Image Ghosting and customer complaints. This procedure needs to be done EVERY TIME software has been reloaded. Line Rate is NOT captured on the SavInfo MOD or in the Health Page Status of the system.

Note:

UNIX/IRIX is case sensitive. Use the commands exactly as shown as "root owner".

1. From the Desktop click on the **C-shell tool**
2. Switch to root owner. Type: **su root** [enter]; (Enter password)
3. Type: **chkconfig fixed72hz off** [enter]

Note:

If you get the following error you may have miss-typed the command. Re-type the command. If it appears the second time, IGNORE IT. Proceed with the step 4.

chkconfig: cannot set flag "fixed72hz": No such file or directory
use "chkconfig -f fixed72hz off" to create and set the flag

4. Type: **/usr/gfx/setmon -x -w 1280x1024_60**
5. Type: **exit** [enter] to exit root
6. Type: **exit** [enter] to exit the C-shell
7. Select **System Shutdown**, (OK to shutdown) (Reboot is essential)
8. Select **Restart**
9. **Choose Signa icon**, (Enter password)
10. Proceed to Section 5 : Signa Host Gamma Setup – NEC 1880SX

5- SIGNA HOST GAMMA SETUP- NEC 1880SX

The GAMMA value is modified to optimize the contrast level of the image mid-tones to more closely represent the same contrast that is filmed.

Note:

At time of System Installation:

- As of July 28, 2003 all NEW LCD monitors will have a Gamma Software CDROM disk INCLUDED in the box, (GE Part # 2376768) which should be used for software versions based on 9.0 and below. Use this CDROM disk and the instructions it contains to setup gamma for those software versions.


Upon replacement of a failed monitor:

Make sure you also order:

- **GE Part # 2376768.** This CDROM disk is labeled "MR Signa® Host Monitor Gamma Correction Software" contains the most up-to date gamma files needed to setup gamma. The software on the CDROM is required to install a replacement monitor. Follow the installation instructions included with the CDROM disk.

WARNING!

MR SIGNA® HOST MONITOR GAMMA CORRECTION SOFTWARE (GE PART # 2376768) SUPERCEDES ALL OTHER GAMMA FILES THAT MAY BE ON VERSION 14 AND 15 LX/8X SERVICE METHODS CDROM.

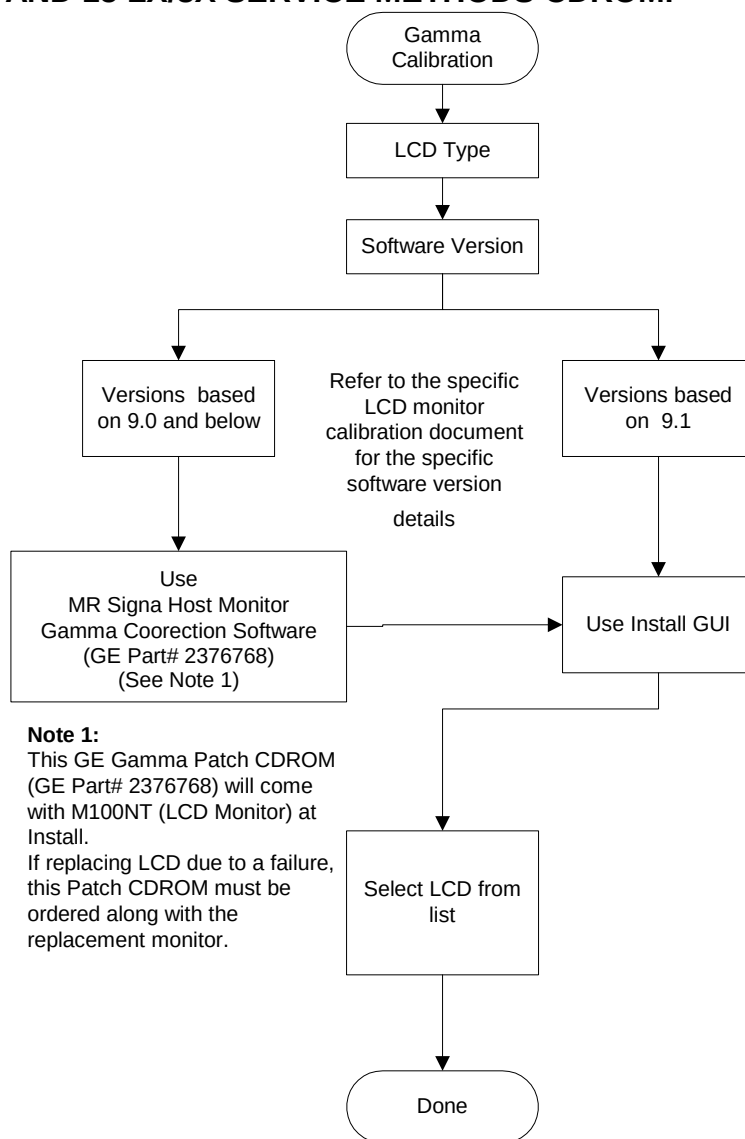

PROCESS FLOW CHART
ILLUSTRATION 1-1

5-1 Gamma Configuration for all SGI Computers

Select the section to use based on your Computer type and Software version, using the table below:

GAMMA PROCEDURES BY SECTION

TABLE 5-1

Computer Type:	Software Version:	Gamma Setting:
All SGI computers	ALL Software versions 9.0 and below. Setting Gamma <u>with</u> the CDROM disk that came with the monitor.	See Section 5-1-1
Octane 1 and Octane 2	9.1, HFO2, HF03 Setting Gamma on monitors <u>without</u> a CDROM Disk	See Section 5-1-2

5-1-1 Setting Gamma with the Patch CDROM – All Software versions 9.0 and below

- If you are Installing the NEC 1880SX monitor for the first time a CDROM labeled: MR Signa[®] Host Monitor Gamma Correction Software is included in the Monitor box. The installation instructions are included with the CDROM disk and should be followed to load the gamma software on the system.
- If you are replacing a failed LCD Monitor, then a CDROM labeled: MR Signa[®] Host Monitor Gamma Correction Software, (GE Part # 2376768) should have been ordered along with the replacement monitor. The installation instructions are included with the CDROM disk and should be followed to load the gamma software on the system.
- This latest CDROM disk should be retained and re-used any time a reload of software occurs on the system.
- If you do not know the software version, type “getver” in a c-shell.


This option is used to configure an NEC 1880SX LCD. To set gamma:

- **If you are replacing the NEC 1880SX with another 1880SX and have not reloaded software:** No gamma changes are necessary. Proceed to Subsection 6 in this document.
- **If you are replacing the NEC 1880SX with the another 1880SX and have reloaded software:** Proceed to STEP 1 below:
 1. Insert the “MR Signa[®] Host Monitor Gamma Correction” CDROM into the Host CDROM/DVDROM drive.
 2. Follow the instructions included with the CDROM to set the gamma for your monitor type.
 3. Proceed to section 6 in this document.

5-1-2 Setting Gamma - 9.1, HF02 and HFO3

This option is used to configure an NEC 1880SXLCD. To set gamma:

- **If you are replacing the NEC 1880SX with another 1880SX and have not reloaded software:** No gamma changes are necessary. Proceed to Subsection 6 in this document.
 - **If you are replacing the NEC 1880SX with the another 1880SX and have reloaded software:** Proceed to STEP 1 below:
1. Open Guided Install and select **System Configure** tab. Refer to illustration 5-1.
 2. Look for the “**Display Type**” field. (See Illustration 5-1)

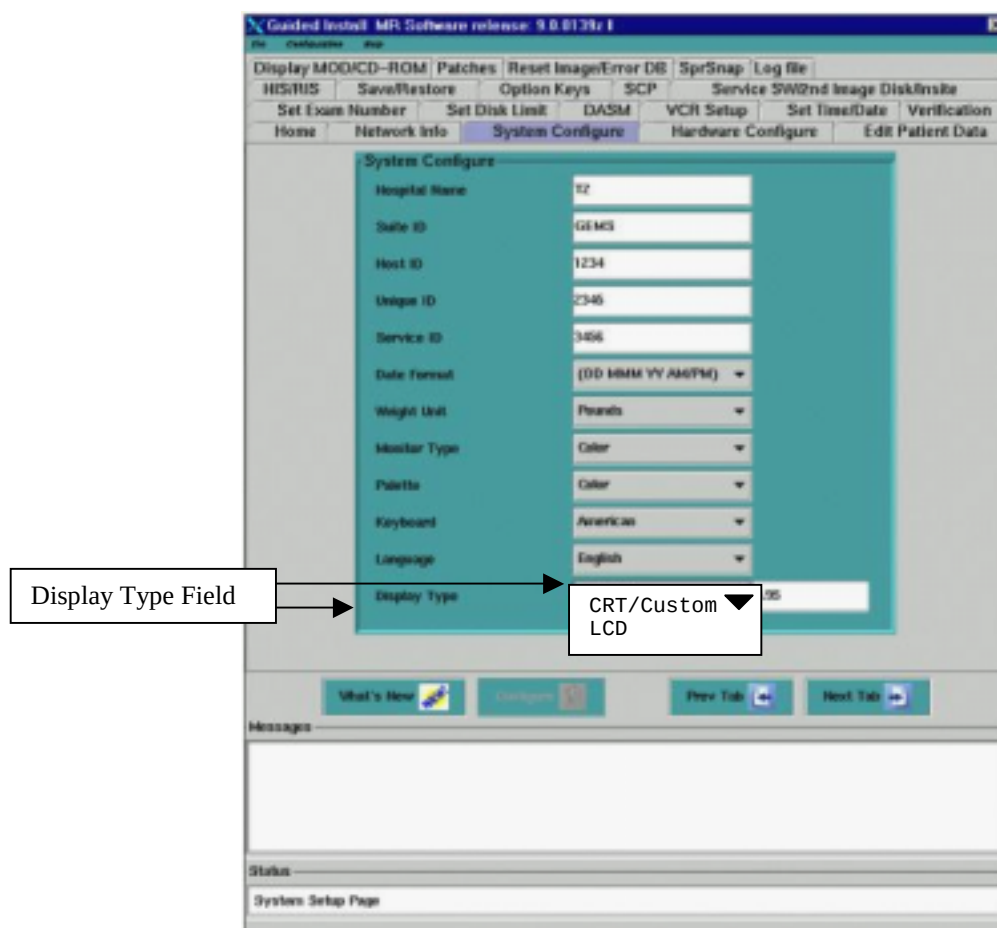

GUIDED INSTALL GUI “SYSTEM CONFIGURE” TAB- DISPLAY TYPE FIELD
ILLUSTRATION 5-1

Note:


The “Display Type” Field has two choices:

- **LCD** – This should be used to configure Gamma for a specific LCD listed in the dropdown box. (Preferred method)
- **CRT/Custom** – This should be used to set a Custom CRT Gamma setting.

Note:



It may be necessary to “Toggle” between the “**CRT/Custom**” and “**LCD**” options in the Display Type Field to select “**LCD**”.

4. If “**CRT/Custom**” is displayed in the Field:
 - a. Click on the “Display Type” Field ▼.
 - b. Drag your mouse over **LCD**
 - c. A Popup Display message will appear, Select [OK] (See illustration 5-2)

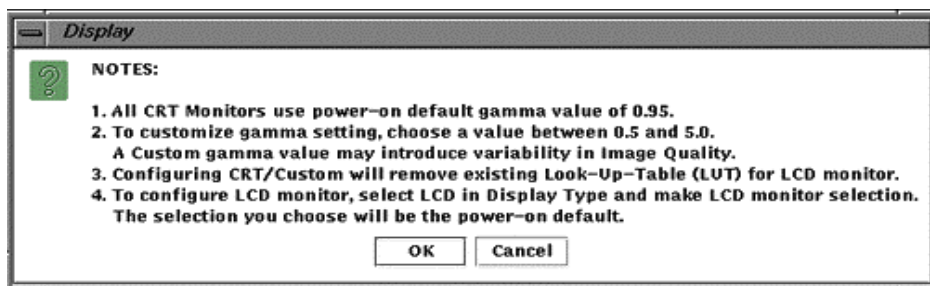

MESSAGE BOX
ILLUSTRATION 5-2

- d. The Field should change to [LCD].
- e. Select **LCD** again. The following menu will appear:

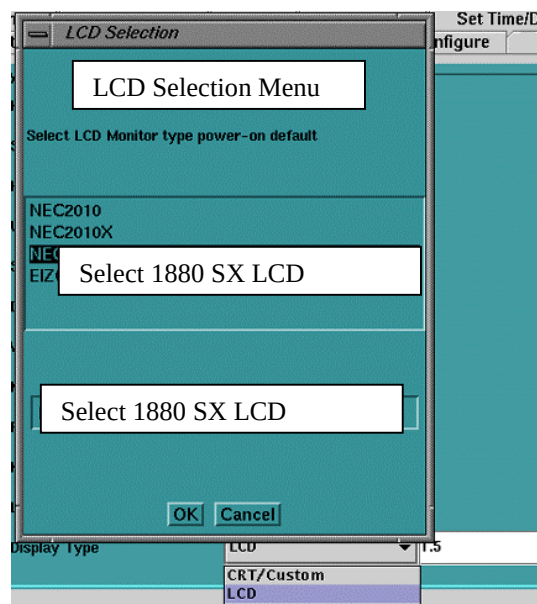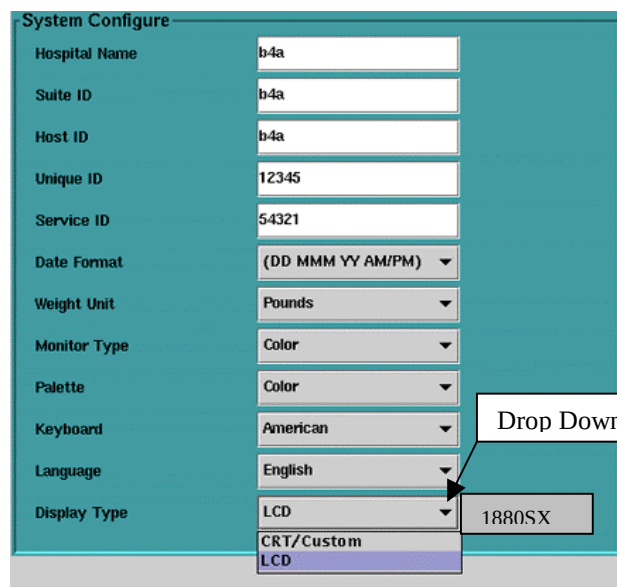


MONITOR CHOICES
ILLUSTRATION 5-3

- f. **Scroll to the correct Display option (1880SX) and release the mouse button.** The display selected will appear in the box to the right.
- g. Click on the [**Configure**] Button on the Install GUI to save the changes. You may see a slight change in the screen intensity when “Configure” is complete.
- h. A Popup message box will appear, Select **OK**.
- i. Select **File** and **Quit** and **OK** from the Install/Reconfigure GUI Window.

- j. **Reboot the System.** A reboot is required to activate the changes. This completes the Gamma Calibration process.
- k. Proceed to Section 6- LCD Auto Adjustment and OSM™ On Screen Manager Operations

6- LCD AUTO ADJUSTMENT AND OSM™ ON SCREEN MANAGER OPERATIONS

This section should get you familiar with the basic operation of the monitor controls. The intent of this section is to adjust the monitor to a "default" condition to proceed with calibrations. It also will allow you to become familiar with the button panels and the OnScreenManager setup.

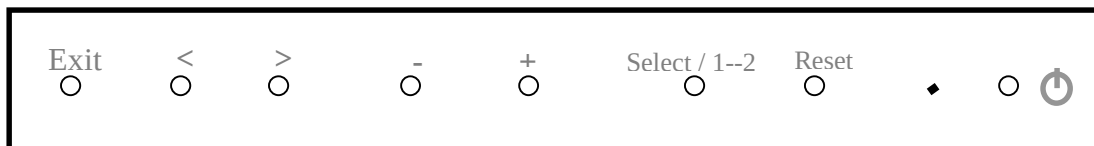

NEC 1880SX LCD FRONT PANEL OSM ACCESS BUTTONS
ILLUSTRATION 6-1


TABLE 6-1
OSM ACCESS BUTTON FUNCTIONS

Exit	Exits the OSM controls, Exits the Main Menu.
< >	Moves Highlighted area left and right within the sub-menus.
- +	Moves bar left/right to increase or decrease adjustment.
Select / 1--2	Enter OSM controls and sub-menus. Activates auto adjust function.
Reset	Resets the highlighted menu to the factory setting.
	On Off switch

6-1 Overview

In this section the following adjustments must be performed:

- 1. Auto Adjust →
- a. Image Adjust (automatic) →
- b. Position Control (Automatic) →
- 2. AccuColor™ set-up (9300k) →

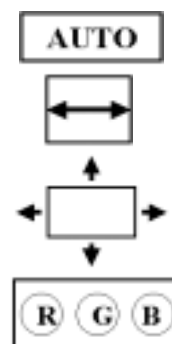


3. Readjust contrast and brightness to customer's recommendation.

6-2 Using the OnScreen Manager

Adjustment bar extended to 80%

NEC 1880SX LCD FRONT PANEL OSM ACCESS BUTTON AND MENU ILLUSTRATION 6-2

1. To access On Screen Manager OSM™, press one of these buttons: +, -, <, >

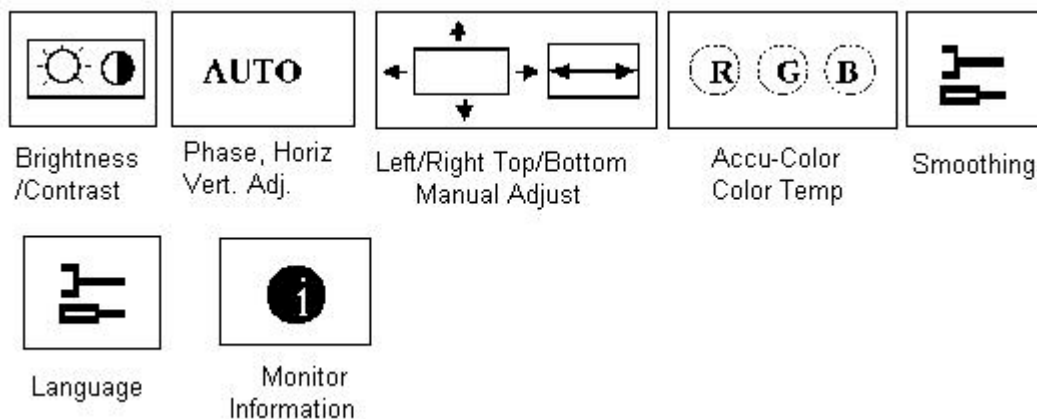

MONITOR ADJUSTMENT ICONS ILLUSTRATION 6-2

Note

To Rotate the OSM Menu Between Landscape and Portrait mode press the OSM button On, then Exit Button, then press OSM button again.

2. Select the **Select / 1--2** Button on the front of the monitor.
2.  Push either the left and right Control buttons - + to initiate the Auto Adjust of the Phase/ Vertical and Horizontal Centering.

Note

Manual Adjustment of the Position Control and Image Adjust H. Size/Fine controls is rarely needed. Automatic adjustment takes care of this.

3. Select the **Select / 1--2** Button on the front of the monitor



4. From the **AccuColor® Control System** icon and use the < > Control buttons to select the first entry row under this icon.
5. Select: **1 (9300)**. (Monitor Color Temperature), Select the **Select / 1--2** Button on the front of the monitor.
6. Leave **Smoothing** setting at **Default**, Select the **Select / 1--2** Button on the front of the monitor.
7. Leave **Language** setting at **Default**. Select the **Select / 1--2** Button on the front of the monitor.
8. Leave **Display mode** setting at **Default**. Select the **Select / 1--2** Button on the front of the monitor.
9. Press the Exit Button. All adjustments are now saved. These settings will remain in effect even if power is removed, or system software is re-booted. Unless manually changed by the operator in the OSM , On Screen Manager.
10. The monitor adjustment is complete. The LCD Display Monitor is now ready for integration to Signa.

Note

For further information, of the OnScreenManager and button identification, refer to the NEC LCD 1880SX User's Manual that came with the LCD Display.


7- FINAL ADJUSTMENT OF CONTRAST AND BRIGHTNESS

The SMPTE (Society of Motion Picture and Television Engineers) test pattern is used to provide a standard image for calibrating the Window and Level of the display for an MR scanner.

7-1 Displaying the SMPTE Pattern

The Window and Level adjustments must now be set to ensure site to site uniformity in image appearance. The SMPTE test pattern is available on the Operator Workstation Host Computer after IRIX, (The operating system) has been booted, and after you have logged into the system.

1. Install and display the SMPTE pattern. Using the steps in illustration 7-1 Below :

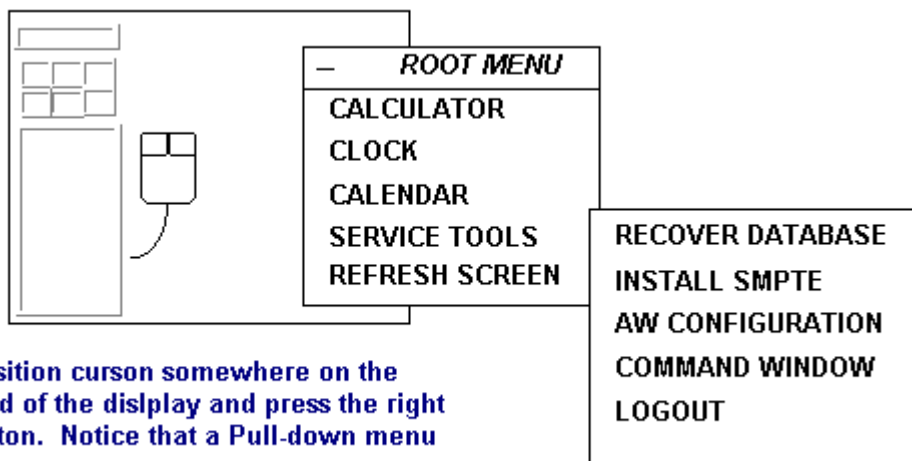

Step 1. Position cursor somewhere on the background of the display and press the right mouse button. Notice that a Pull-down menu appears.

Step 2. Slide the mouse down to [Service Tools], and then over to [Install SMPTE]. This will load the SMPTE pattern into the database as Image 1000.


INSTALLING THE SMPTE PATTERN ILLUSTRATION 7-1

2. Displaying the SMPTE Pattern by following the steps in illustration 7-2 on the next page. Switch to the Browser and select the SMPTE pattern from the patient list and display it.
3. If necessary, use the FORMAT option to select Full Screen mode.

Step 1. Point to and click on the Display icon. Notice that the Browser comes up.

Step 2. After the Browser comes up, use the scroll bar on the right side of the display to find Image 1000 SMPTE.

Step 3. Point to and single-click on the SMPTE entry.

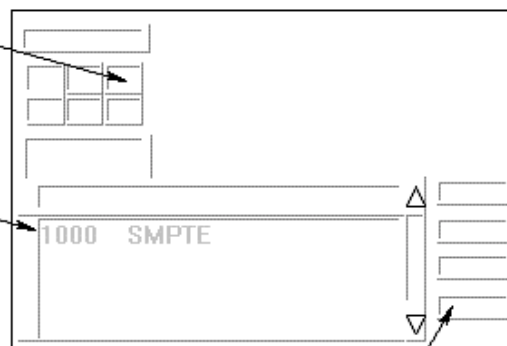


Step 4. Point to and single-click on the Full Viewer to display the SMPTE so that it is full screen size.

Step 5. To get back from displaying the SMPTE pattern, hit <ESC>.

INSTALLING THE SMPTE PATTERN ILLUSTRATION 7-2

7-2 Adjusting for proper Window and Level Setting

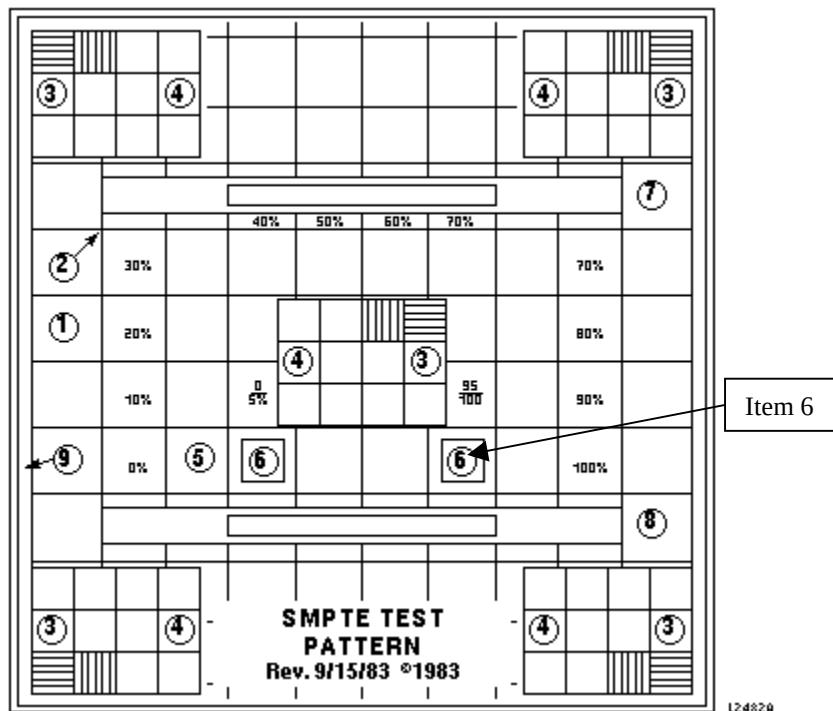

SAMPLE SMPTE TEST PATTERN
ILLUSTRATION 7-3

1. With the cursor inside the displayed image, hold down the middle mouse button and move the mouse in the horizontal plane. View the window control value at the base of the image and **set window control value to 100**.
2. With the cursor inside the displayed image, hold down the middle mouse button and move the mouse in the vertical plane. View the level control value at the base of the image and **set level control value to 1024**.

7-3 Fine Tuning the Contrast and Brightness using the SMPTE Pattern

In this section you will adjust the Contrast and brightness controls of the LCD Monitor to optimize the SMPTE Test Pattern for the Window and Level controls.

Minimal adjustment should be necessary at this time. However, it is important that you "fine tune" the Contrast and Brightness level using an industry Standard Test pattern.

This is an iterative process by adjusting the contrast, then brightness, and then back to contrast.

Note

To reduce visible tearing or smearing of the pattern, LOWER the contrast.

1. To access On Screen Manager OSM™, press one of these buttons: +, -, <, >

2. Use the **Control** buttons - + on the front of the monitor, to adjust the Brightness and Contrast .
3. Look at the Pattern. If the 5% and 95% patches are not visible.(item 6, Illustration 7-3), adjust the brightness control of the display monitor until you can just barely see them. You may also need to re-adjust the contrast, if tearing or smearing of the pattern occurs. (items 1 and 2, and 5 through 8, Illustration 7-3).Continue adjusting both brightness and contrast until the image is as crisp and clear as you can make it, and the 5% and 95% patches are just barely visible.
4. When satisfied with the display, exit the menu to save Contrast/Brightness settings.

Note

These settings will remain in effect even if power is removed or the system is re-booted. Because an operator can tamper with the Brightness and Contrast Adjustments, It is highly recommended that you make all operators aware of the affect they will have on system image quality and the camera images.

Note

If the operator chooses to deviate from this Brightness and Contrast adjustment, ensure the camera imaging and calibration has NOT been affected.

5. Remove the SMPTE Pattern by selecting a different Exam from the browser.

8 - CAMERA CALIBRATION

This procedure describes the steps necessary to verify and set the camera parameters. Once the display is re-calibrated, it is essential to re-calibrate the camera before the system is used for filming. Although a qualified GE Service Engineer could perform the steps below, it is recommended that the following procedure be performed with the on-site assistance of the camera vendor field engineer.

8-1 DASM Interpolation Setup

1. Select the **Install** soft key from the Service Desktop. When prompted, login using the password **operator**.
2. Select the **DASM** folder from the top of the Guided Install GUI once it appears.
3. Set the DASM Interpolation method to **linear**.
4. Before exiting the GUI, insert the SaveInfo MOD or new MOD into the drive. Go to the top left corner of the Install GUI and select **<File>** and then **<Save GI Configuration to MOD>**. This process does not create a SaveInfo disk. It just creates a copy of the information already entered in the GUI tabs for use in the next software install.
5. Exit the install GUI by selecting **<File>** then **<Quit>**. If any error conditions still exist, you will be warned that no changes will be made before exiting.

8-2 Camera Imaging Look-up Table

Verify with the camera vendor field engineer that the currently installed camera lookup table is designed to provide perceivably linear contrast for the light box conditions in the customer's viewing area. If not, request that the camera vendor field engineer replace it with the appropriate look-up table.

8-3 Camera Maximum Optical Density

Note that for optimal reviewing, the light box luminance of the diagnostic region of the film should be in the range of 50 to 500 nits. Vary the maximum optical density setting of the camera to compensate for the light box and to meet this value. A good starting position is a maximum density of 2.8.

Note

The final OD settings may be refined by the radiologists performing the image review.

8-4 Camera Contrast

1. With the maximum /minimum optical densities set to compensate for the review area's light box, select a look-up table for your camera that will produce a perceivably linear gray scale for the same light box and the overall ambient light conditions of the viewing area.

Note

The DICOM 3.14 Standard specifies the Barten's curve for linear perception. It is recommended that the manufacturer base perceptual linearity on this curve.

2. Film the SMPTE pattern on a 1-on-16-format display. Verify that the 5% and the 95% levels are visually equivalent. If not, perform a Contrast test with the SMPTE pattern. Select the new contrast setting from the contrast image set. A good value for the Imation DryView is 3. Use the camera's calibration procedure to set the contrast setting. Ensure that the camera maintains a perceivably linear gray scale.

Note

Filming the SMPTE pattern for contrast calibration may be optional for the camera vendor.

3. Ask the technologist to display a clinical image and set window and level controls for the desired appearance. A good image to start with is a sagittal or axial head image.
4. Capture the image on the keypad or host control interface.
5. Print a Contrast Test film. Ask the technologist to select the image that best matches the displayed image on the monitor.
6. Observe the image number below the selected image and set the Contrast control to this value.

8-5 Anatomical Filming

This portion of the procedure requires the technologist to verify the camera settings with true anatomical images.

1. Film representative anatomical images to confirm the settings. The image set should include T1 and T2 head images, joint images and c-spines.
2. Observe the accuracy of the low-tones, mid-tones and high-tones. If a filmed image is found to not be equivalent then re-calibrate the camera based on the customer's evaluation.

9 - TROUBLESHOOTING GUIDE

The following section provides suggestions for troubleshooting the LCD monitor:

9-1 Explanation of Default Gamma Settings

The Default gamma setting is different for different software versions. The Table below indicates the suggested nominal settings. The Software Default Settings are set at boot up. When a Look-Up-Table is installed, these default settings are ignored. If the "setgamma" command is used, the "colorPix.cfg file and the boot script (Startmon) is updated, with the new single value.

TABLE 9-1
DEFAULT GAMMA SETTINGS

LCD Monitor	Nominal Gamma value
NEC 1880SX	0.79

As the Gamma setting gets lower, the display will become dimmer, (less contrast).

9-2 Complaint: Poor SNR, Dynamic Range Lack of Grayscale.

- It is important to note that the software default gamma setting of 1.5 causes extremely poor image quality on both the LCD and CRT. Customers complain of Ghosting or Poor SNR or lack of dynamic range in the grayscale. Before troubleshooting any system problem related to these symptoms check the Monitor Gamma Calibration, or make sure you have the correct Monitor Selected.
- Gamma Calibration tables (LUT's) are unique from one monitor type to another. DO NOT use a gamma table designed for a different monitor. If the gamma table or Monitor Type does not exist, start with a single gamma value of .95. It may be necessary to "tweak" this setting for your customer. Use the "setgamma" command to do so.

9-3 Complaint: Ghosting or Blurring images or Sheering alphanumeric LCD 's ONLY.

- Ghosting and blurring related to areas of gray in an image. This can be related to system problems but also can be attributed to monitor Frequency Rate (Frame and Swap rates). ALL LCD's should be set to 60Hz.

Checking the Frame Rate and Swap Rate Setting:

1. From the Desktop click on the **C-shell button**
2. Switch to root owner. Type: **su root** [enter]; (Enter password)
3. Type: **toolchest&** [enter]
4. Click on **Find → Control Panels**
5. **Double click on Display Icon** to open Display tools.
6. Verify: **Active Size (Pixels): = 1280x1024 and FRAME RATE and SWAP RATE =60.00 Hz** (Refer to Illustration below:



FRAME RATE AND SWAP RATE
ILLUSTRATION 9-1

7. If it is the Frame Rate and Swap Rate are not set to 60Hz:
8. Select **1280X1024 60** from the left menu.
9. **Click on "Load" button.**
10. **Click on OK** when the Pop-up box asks you if you want to Load the setting now.
11. **Click on OK** when the pop-up box asks you if you want to keep the Format.
12. **Click on OK** when the pop-up box asks you if you want to make this setting the Power-On Default.
13. Close the Graphics Backend Window: Click on **File → Exit**
14. Close the Control Panel Window: Click on **Page → Exit**
15. Close Toolchest Window: **Right Click and hold on the word Toolchest** at the top of the toolchest window.
16. Drag the mouse to **Close**
17. In the C-shell type: **exit** [enter] to exit root
18. Type: **exit** [enter] to exit the C-shell.
19. Reboot the system.

9-4 Camera Film has changed after Monitor replacement

Note:

If there is no remote camera in the list, you will not be able access the Film Composer Menu.

- Make sure that the Monitor and Camera Calibration procedures have been performed completely and the correct Monitor Type has been selected and installed.
- Then try setting the camera to "**Replicate**" instead of linear mode, or vice versa.
- For Kodak Cameras try setting to "**Bi-linear**" mode.

- For Other cameras try setting “**Replicate**” mode.
Procedure:
 1. Open the Image Browser
 2. Open Film Composer
 3. Select **configure**
 4. Select remote printer being used by the site.
 5. Select **Update**
 6. Under magnification Type:
Select “**Replicate**” or “**Bi-Linear**” option. (Kodak try “Bi-linear”)
Select **OK**
 7. Select **Done**
 8. Select **Quit**
 9. Select **OK**
 10. Return to original desktop.
 11. Check film and adjust camera if necessary.

9-5 No Picture

- The signal cable should be completely connected to the display card/computer.
- The display card should be completely seated in its slot.
- Power Switch and computer power switch should be in the ON position.
- Check the signal cable connector for bent or pushed-in pins.

9-6 Image Persistence

- Image persistence is when a “ghost” of an image remains on the screen even after the monitor has been turned off. Unlike CRT monitors, the LCD monitor’s image persistence is not permanent.

Note

It is still recommended that a screen saver be used whenever the screen is idle.

9-7 Image is unstable, unfocused, or swimming is apparent

- Signal cable should be completely attached to the computer.
- Check the monitor and your display card with respect to signal timings.
- (LCD’s Only) Change the video mode to non-interlace and if possible use a 50 Hz refresh rate.
- Adjust Monitor controls to focus and adjust display. Refer to the vendor manual that came with the LCD Monitor for the operation of this function
- (LCD’s Only) If your text is garbled, change the video mode to non-interlace and if possible use a 50 Hz refresh rate.

9-8 LED on monitor is not lit (no green or amber color can be seen)

- Power Switch should be in the ON position and power cord should be connected.
- Make certain the computer is not in a power saving mode.

9-9 Display image is not sized properly

- Push the Auto adjustment button.

9-10 Selected resolution is not displayed properly

- Push the Auto adjustment button.

9-11 Diagnostic Image Quality has been traced to problem with LCD setup.

- Operator has adjusted the Brightness and Contrast of the monitor which has affected Camera/Film imaging. Re-calibration may be necessary. Start with Section 1 of this procedure.

9-12 Display image has a green cast to it

- Select “**TYPE**” in the OSM™ Information menu and press the < or > control button.
- Set the Color Temperature to 9300K

9-13 Display image ghosts or gray lines in areas of extreme black & White transitions.

- The LCD Refresh Rate is too fast and must be slowed down to allow for the Pixels to react properly. To work around the issue all LCD monitors should be set to Refresh Rate, (Also known as a Line Rate or Frame and Swap Rate) of 60Hz. . The trade-off will be: Better Black & White Images at the sacrifice of color smoothness. Letters such as “V” and “W” will appear more jagged. See Section 4 of this procedure.

9-14 Display image Blurry or Blotchy images Described as poor grayscale.

LCD's are highly dependent on gamma correction. Perform the gamma correction procedure to solve this problem. See Section 5 of this procedure.

Note: The off-the-shelf LCD's on Signa are optimized for PC use from the factory, NOT diagnostic Imaging. The LCD Monitors are NOT “Plug and Play” with regards to the Signa Applications Software and SGI Computer Video system. This requires that a complete monitor Setup and calibration be performed at:

- Installation
- During a replacement of a failed monitor.
- Any time Software is reloaded.

9-15 Customer complaints of poor grayscale, Window or level, lack of definition.

LCD technology uses three color film plates and a backlight to display different shades of color and gray. They are highly dependent on gamma correction. The off-the-shelf LCD's on Signa are optimized for PC use from the factory, NOT diagnostic Imaging. Attaching an LCD to the Signa system without calibrating it will cause image quality problems that can be interpreted as system problems, leading to unnecessary:

- Hardware replacement
- Protocol Troubleshooting
- Customer Satisfaction Issues.

It is important to note that LCD Monitors are NOT “Plug and Play” with regards to the Signa Applications Software and SGI Computer video system. This requires that a complete monitor Setup and calibration be performed at:


- Installation
- During a replacement of a failed monitor.
 1. Any time Software is reloaded

9-16 LCD FAQ's:

- The Monitor calibration information is NOT saved during the SavInfo backup routine or accessible via Healthpage.
- There is a big difference between CRT and LCD monitors. CRT's are much less sensitive to a wrong gamma setting.
- Gamma Tables (Look-Up-Tables or "LUT's") are highly specific to the LCD model and should only be applied to the correct model they were designed for.
- Gamma Tables consist of thousands of luminescent data points that are used by the electronic circuitry to generate uniform transitions from one color or grayscale level to the next.
- All recent LCD's used in Signa have a "best fit - Single gamma value" as well as a complete Gamma Table to maximize flexibility when setting up the monitor.
- It is recommendd that the Gamma Table always be used as a default.

9-17 What is Gamma?

Simply stated gamma normalizes the grayscale so each level is unifor in size and shape across the entire display:

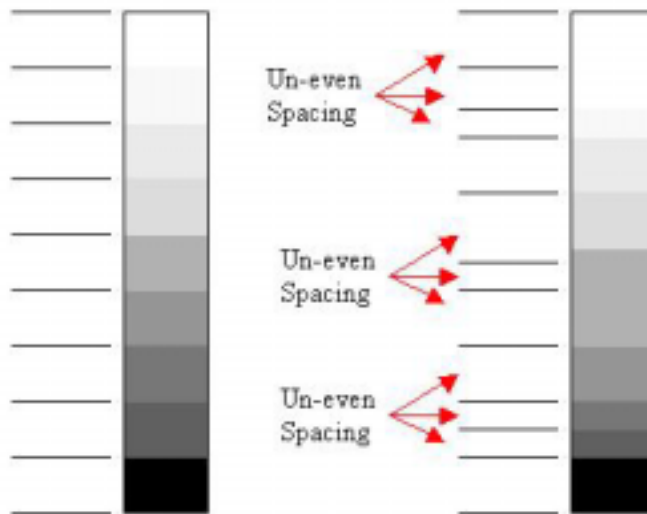

GRAYSCALE COMPARISON
ILLUSTRATION 9-1

Without getting into complicated equations "Gamma correction" adjust the grayscale or color gradient so each "step" is exactly the same size and shade of color or gray. If the "steps" are not uniform in size and color, the monitor will display the error as poor transitions from one step to the other. This appears to the eye as small ghosts or ragged edges.

9-18 How do I set the correct gamma?

- When the Monitor Type is selected in the Install GUI the correct gamma table (LUT) for that monitor is automatically loaded. The gamma LUT remains active even if a reboot occurs.

9-19 Can I use the gamma command to set gamma?

- NO. This command is no longer an effective way to change gamma permanently. This command should only be used if the need arises to temporarily change to a single gamma value, Overriding any previously installed Gamma LUT.
- In the past, we used to adjust gamma using this command. THIS HAS CHANGED. Using the gamma command on systems with an LCD will most likely CAUSE IMAGE QUALITY PROBLEMS.
- Using a Gamma Table will in most cases produce better display image quality than a single gamma value. A Single gamma value is much less effective on LCD's at providing a uniform grayscale.

9-20 How do I determine if the correct gamma Table or setting is installed?

- Open up the Install GUI. Look at the Display Type Field. To the right of the Option box it will either display a monitor type or a single number.
 - o If a Monitor Type is in this box then a "Gamma Look-Up-Table" for that monitor is installed.
 - o If a "Single Value" is displayed, then the CRT/Custom Option was used to set a Single gamma value.

9-21 When I type "gamma" in a c-shell I see a number. Is that the gamma setting?

- NO. It is simply a "Meaningless number" that no longer holds any value. When a gamma table is installed, (By simply selecting the monitor Type in the Install GUI) This value is meaningless. It is a left over from earlier versions of software.

9-22 When I type "gamma 1.05" in a c-shell, the LCD changes in appearance. Why?

- This is complicated because IRIX Operating System Commands such as "gamma" take priority over the Signa Applications Software. The Install GUI Gamma Tools are written in IRIX Script and have less priority. Since the O/S gamma command was invoked with a value, that value will now OVERRIDE the Applications software. The Single Value of 1.05 temporarily overrides the gamma LUT until a reboot occurs. (Upon reboot, the Gamma LUT resumes control).

9-23 A gamma value of 1.05 is what we always used in the past can we still use it?

- The gamma value of 1.05 is no longer used. All LCD and CRT Monitors now use .95 which was determined to be the best default "Single value".
- DO NOT USE THIS COMMAND ANY MORE. The process HAS CHANGED. Using the gamma command on systems with an LCD will most likely CAUSE IMAGE QUALITY PROBLEMS and filming problems.

- Using a Gamma Table rather than a single gamma value, (Such as 1.05) will in most cases produce much less IQ problems than a single gamma value.
- A Single gamma value is much less effective on LCD's ,at providing a uniform grayscale.

9-24 My customer does not like the gamma setting that the Gamma Table provides.

- By all means change it! use the Install GUI, under Display Type" Select "CRT/Custom" and type in a "Single Gamma Value" you and the Customer agree on is satisfactory. A reboot will make the Single setting permanent.
- There is a Single optimized "Nominal" Gamma Value for each monitor type as well as a gamma table. Consult the Troubleshooting sections within the proper monitor calibration procedures to find the "Nominal" value for the monitor type.
- Note: Using a Single value affects the camera to a much higher degree than a Gamma Table does, so be prepared to adjust the camera as well.

9-25 Single value optimization technique.

- To minimize the time needed to find a Single gamma value that is acceptable by your customer, you can still temporarily use the IRIX O/S gamma command in a c-Shell to "Tweak" a single gamma value.
 - o Typing "gamma *n.nn*" (where *n.nn* is the nominal single gamma value for the monitor type. See 9-22) It will generate an immediate change to the monitor. (Increasing the value will increase brightness and contrast together. Decreasing it will have the opposite affect).
 - o Continue to adjust the gamma setting, checking images and films until you arrive at a setting you like. Note the final value.
 - o DON'T FORGET THIS STEP: To make the value permanent Open the Install GUI and under Display Type" Select "CRT/Custom" and type in the new "Single Gamma Value" .
 - o Reboot the system and this will now be the permanent setting.

Note: If you open a c-Shell and type "gamma" it may not match what you installed in the Installation GUI. This is normal. See 9-22).

A- APPENDIX A

A-1 Focus, Screen position, Clock Rate and Phase Adjustments.


The Auto Adjust procedure done in Section 5 takes care of these for you. All other values in the Screen Manager should stay at the default values except Contrast and Display.

Use this appendix only if you see problems with Focus, screen position relative to monitor mask or phase problems (horizontal bars sweeping through the LCD). Or any vertical problems.

A-2 Adjusting Focus Using LX Calibration Tools From Toolchest (Alternate)

This section is normally not needed. However, if you cannot adjust the monitor satisfactorily, with the front panel buttons and OnScreenManager, you may want to try the adjustments described in this section.

1. From the Tools Menu on Signa Open a **C-Shell**.
2. Switch to root owner. (Type: **su** <Enter> password: **operator** <Enter>)
3. Type: **[toolchest&]**. A pulldown menu will appear in the upper left-hand portion on the display.

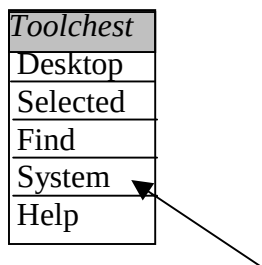

TOOLCHEST MENU
ILLUSTRATION A-1

4. Select **[System]** with the Right mouse button.
5. Select **[Confidence Tests]** from the extended menu. A "confidence test" pop-up will appear after a few seconds.
6. From the "confidence test" pop-up, select the **[Monitor Icon]**. A monitor calibration menu will appear on the screen.
7. Select **[Focus]** from the "Test Options" menu by using the **Tab** key.

8. Put the mouse pointer on the "Test Options" Menu. Click the **Middle mouse button** to toggle (Hide) the "Test Options" menu.
9. Push the auto-adjustment button on the front of the monitor. This should setup Focus, Horizontal, Vertical Centering, sync to the proper Clock Rate and monitor Phase.
10. Observe the test pattern displayed, it should now be centered and in focus. This should be all that is necessary to adjust the display controllers' output to match the LCD monitor.
11. Toggle the "monitor calibration" toolbar back on by selecting the middle mouse. Select the grayscale menu option.
12. Toggle the menu back off by selecting the middle mouse again.

A-3 Setting up LCD Display for Contrast/Brightness. Using LX Calibration Tools From Toolchest (Alternate)

This alternative procedure uses the Monitor Calibration tools and patterns. Steps 1 through 4 should be done while observing the image for crisp lines, overall uniformity of pattern, and to setup the LCD monitor so it is at its brightest level without "Blooming the image or adding distortions". The intent of this part of the appendix is to achieve the best image appearance you can with the LCD. At the completion of this section you should check the display against industry standard SMPTE pattern. (See Section 7 of this procedure - LCD CALIBRATION FOR OPTIMUM IMAGE VIEWING)

1. Select "**Contrast/Brightness**" button and **push enter**.
2. Adjust brightness to maximum level **100%**.
3. Look at the monitor and adjust the Contrast and Brightness to optimize the test clarity of the test pattern. This is an iterative process between both settings.
4. When satisfied with the displays' contrast and brightness setting is properly achieved, **push the enter button to save the adjusted Contrast/Brightness settings**.
5. Exit the Monitor Calibration by selecting "**Quit**".
6. Select "**File**" on the Confidence Test Window and choose "**Exit**".
7. **Click on the word "Toolchest"** at the top of the menu with the **Right mouse button** and select "**Close**".

8. Exit the C-Shell by typing: **Exit [Enter]**

A-4 Monitor Adjustments

The monitor should be positioned no closer than 16 inches and no further away then 28 inches from your eyes. The optimal distance is 24 inches for either of the monitors.

Note


Allow the monitor to warm-up for 20 minutes before performing any adjustments.

Note

The response time for the LCD monitor is 100 ms. It is normal to see a “trail” on the screen if the mouse is moved quickly.

A-5 Contrast and Brightness Adjustment For NEC1880SX LCD Monitor

You will need to use the SMPTE pattern for this adjustment. Refer to that section in this procedure.

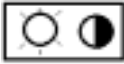
1. From the Contrast and Brightness icon  use the **Control** buttons to adjust the contrast lower if the maximum setting causes visible tearing or smearing of the pattern or alphanumeric characters (items 1,2,5-8 in Illustration 7-3).
2. From the Contrast and Brightness icon  use the **Control** buttons to adjust the brightness of the display monitor if the 5% and 95% patches are not visible (item 6 in Illustration 7-3).

Note

It may be necessary to re-adjust the contrast if tearing or smearing of the pattern or alphanumeric characters occurs (items 1,2,5-8 in Illustration 7-3).

A-6 Control Descriptions- NEC 1880SX LCD Monitor

The NEC 1880SX LCD Panel Display does not have an RGB (BNC) style input. The video cabling uses an S-VGA type connection. It also has additional DVI/D-SUB input capability. This feature allows for the switching of inputs between two LCD Monitors, (Not Used by GE).

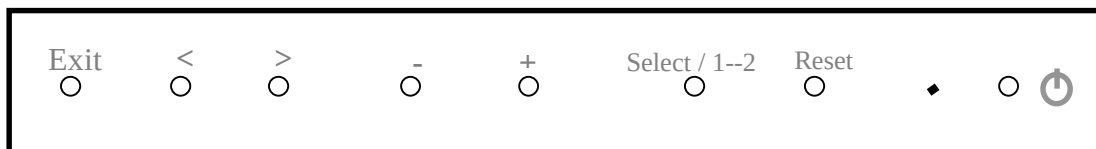

BUTTONS ON FRONT PANEL
ILLUSTRATION A-2

Note

Adjustment information is not lost when the power is removed from the monitor.

A-7 Front Panel Controls

The functions of the OSM™ controls on the front of the monitor are described in Table A-1. To access the OSM, press any of the control buttons (+, -, ◀, ▶) or the RESET/OSM or EXIT buttons.

Note

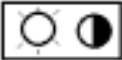
When RESET is pressed, a warning window will appear allowing you to cancel the reset function.

TABLE A-1
FRONT PANEL CONTROL FUNCTIONS - MODEL 1880SX LCD MONITOR

Exit	Exits the OSM controls, Exits the Main Menu.
< >	Moves Highlighted area left and right within the sub-menus.
- +	Moves bar left/right to increase or decrease adjustment.
Select / 1--2	Enter OSM controls and sub-menus. Activates auto adjust function.
Reset	Resets the highlighted menu to the factory setting.
	On Off switch

A-8 On-Screen Controls

This section describes the adjustments available for the OSM controls. The controls labeled with “•” are required during the set-up of the monitor.

- **Brightness and Contrast** → 

Brightness: Adjusts the overall image and background screen brightness.


Contrast: Adjusts the image brightness in relation to the background.

Auto Adjust Contrast: (Analog Input Only) Adjusts the image displayed for non-standard video inputs.

Auto AUTO BRIGHTNESS (Analog input only)
BRT

- **Auto Adjust** → 

Automatically adjust the Image Position or the H. Size setting.

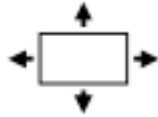

Position Controls →

H. Position: Controls Horizontal Image Position within the display area of the LCD.

V. Position: Controls Vertical Image Position within the display area of the LCD.

Auto: Automatically sets the horizontal and vertical image position within the display area of the LCD.



• **Image Adjust Controls →**

H. Size: Adjusts the horizontal size by increasing or decreasing this setting.

Fine: Improves focus, clarity, and image stability by increasing or decreasing the Fine setting.

Auto This feature automatically adjusts the fine setting for changing signal conditions.


FINE (Analog input only).



• **AccuColor® Control System →**

Six color presets select the desired color setting. This setting must be set to selection 1 - **9300**. If a setting is adjusted, the name of the setting will change to Custom. NATIVE preset cannot be changed. It is the Default setting for the LCD.

R, B, G: Increases or decreases the color depending upon which is selected. The change in color will appear on screen and the increase or decrease will be shown by the color bars.

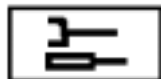

Tools 1 →

You can choose where the OSM control image appears on the screen. Selecting OSM Location allows you to manually adjust the position. Normally, no adjustment is necessary. Use default values.

a-a Sharpness- Adjustment to soften the image.


Smoothing - (Normal) Default

Expansion Mode- (Full Screen) Default

Video Detect- (First Detect) Default

DVI Selection- (Default) No effect on S-VGA systems.

OFF-TIMER- Monitor turns off after predetermined setting.

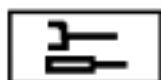

Tools 2 →

You can choose where the OSM control image appears on the screen. Selecting OSM Location allows you to manually adjust the position. Normally, no adjustment is necessary. Use default values.

Language- English

OSM Position- Default location

OSM Turn OFF- Default

OSM Lock Out- This control completely locks out access to all OSM control functions. When attempting to activate controls while in Lock Out mode, a message screen will appear indicating that the controls are locked out. Lockout is NOT recommended.

OSM Rotation- Rotate OSM from landscape to Portrait.

Resolution Notifier- Optimal Resolution is 1280x1024. If ON is selected a message will appear on the screen 30 seconds, notifying the user that the resolution is NOT at 1280x1024.

Xyi - Resolution Notifier- This optimal resolution is 1280x 1024. If ON is selected a message will appear on the screen 30 seconds, notifying you that the resolution is NOT at 1280x1024.

Hot Key- You can adjust the Contrasts and brightness directly.

Factory Preset- This allows you to reset all OSM control settings back to the factory settings. Highlighting the setting and then pressing the RESET button can reset individual settings.

Information → 

Indicates the current display resolution, frequency setting, and type of Sync signal of the monitor. All items under this menu should not be changed. Use default values.

Note

Mode Change should only be used if the monitor does not recognize a resolution. This shouldn't ever be the case with our system. You can change to the appropriate resolution by selecting the Mode information and selecting (increase or decrease) the corresponding option.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	Sept 24, 2002	D. Hofstetter	Initial Version.
1	Dec 2, 2002	D. Hofstetter	Added information to Section 9-Troubleshooting
2	Feb 19, 2003	D. Hofstetter	Added Indigo2 configuration and further clarification about the gamma settings and LUT's.
3	Feb 26, 2003	D. Hofstetter	Added Troubleshooting tips and clarified procedure based on Software version.
4	May 12, 2003	D. Hofstetter	Added Verification section to Section 4 –Setting the Frequency Rate
5	June 4, 2003	D. Hofstetter	Changed Line Rate verification procedure.
6	July 9, 2003	D. Hofstetter	Added section that refers to setup using the Gamma Patch CDROM.
7	July 25, 2003	D. Hofstetter	Removed info regarding old gamma patches and added Gamma Troubleshooting info..