

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1- HOST RECONFIGURATION OVERVIEW	2
1-1 New Installations and Upgrades to Signa LX	2
1-2 Recovery for system with YP(NIS) Service set ON by mistake.....	7
2- MONITOR GAMMA VALUE RESET	8
3- LOAD SERVICE SOFTWARE.....	9
3-1 Load Restricted Service (Only available for GE use.).....	9
3-2 Load Advanced Service (Only for sites with a valid Advanced Service Package Limited License.)	10
5- PERFORM SAVE INFO.....	11
6- FINAL STEPS.....	11
REVISION HISTORY	12

Note

This procedure applies to any system Release 9.0, 8.4 (ACGD), ASP1, and ASP2.

1- HOST RECONFIGURATION OVERVIEW

All Octane Host computers for new installations are pre-loaded with software that is pre-configured for most default values. You must reconfigure the Host computer for the correct configuration files using the Guided Install GUI.


For systems that upgrade hardware, or change the configuration, the Guided Install GUI is also used to change the software to work with the new configuration.

1-1 New Installations and Upgrades to Signa LX

The Guided Install GUI is used to reconfigure the software to match your customer's environment and the hardware that is actually shipped to your customer. For a new installation review data in all tabs and update any information to match actual hardware, system, network ,or other configuration information.

If this is a new/upgraded Signa installation or if the system just needs to be re-configured, do the following to start the Configuration GUI:

1. With Signa up and running, go to the Service Desktop. Click the **[Guided Install [GI:FE Mode]** button then **[START]**. A C-Shell will open and ask for root password ("**operator**" <Enter> is the default). The Install GUI will load and start. See Illustration 1-1.

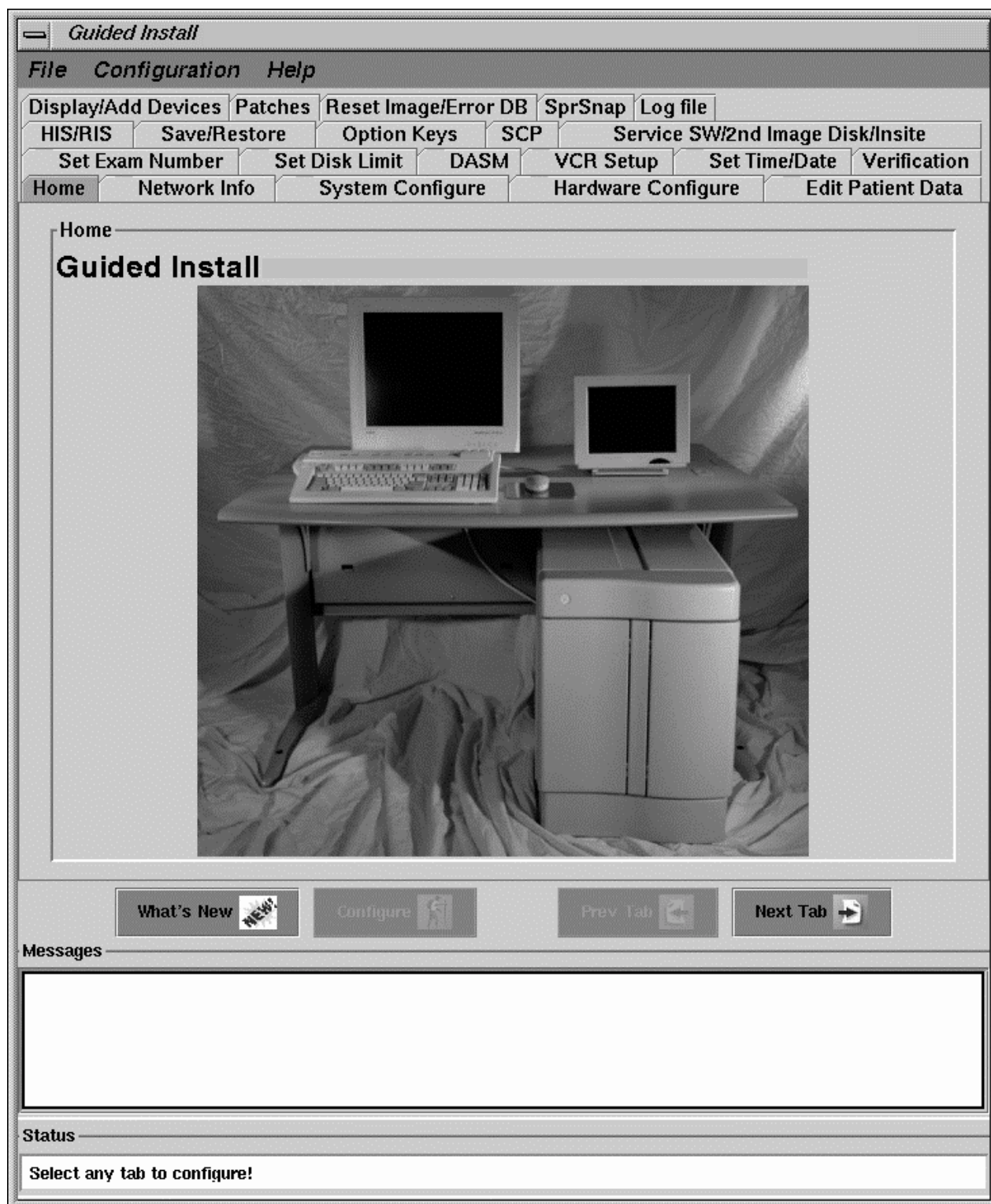

GUIDED INSTALL - HOME SCREEN
ILLUSTRATION 1-1

2. Select **[Next Tab]** button on the Home Screen to bring up the **Network Info** tab. See Illustration 1-2.

GUIDED INSTALL – NETWORK INFO TAB

ILLUSTRATION 1-2


3. Help documentation is built into the GUI. Push the **[Help]** button to view the help information to be used to fill out the required information.

Caution

Leave YP(NIS) Service set to OFF, unless you are 100% sure that there is a hospital network with Yellow Pages service that the customer wants to access. Setting YP(NIS) Service to ON when it should not be can prevent Signa from being able to boot up. See Section 1-2, “Recovery for system with YP(NIS) Service set ON by mistake” if needed.

4. Push the **[Next Tab]** button to move through tabs in order. Modify entries to match your system as required. Continue to use **[Help]** to provide guidance specific to each tab.

5. Note: the **Time/Date** tab is the only tab that the Configure button will never be highlighted to set options. There are separate buttons for configuring this tab (see Illustration 1-3). If you need to re-configure the Host time and/or date at any time, do the following:
- Select the Time Zone your system is located in, then push the **<Set Time Zone>** button.
 - Select the date and local time, then push the **<Set Time/Date>** button.

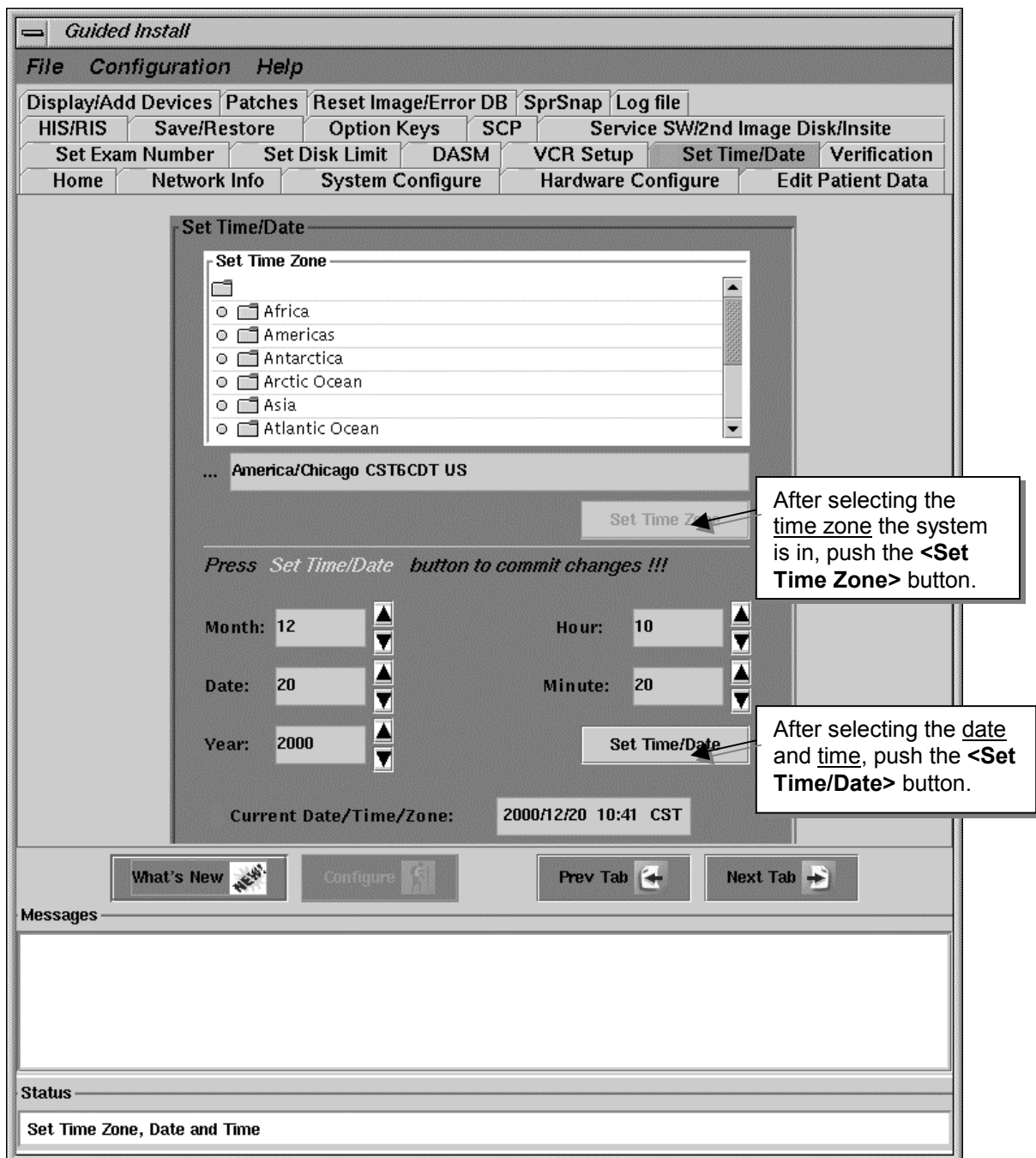

SET TIME/DATE MENU
ILLUSTRATION 1-3

6. There are three regulatory bodies that have defined limits for the use of MR devices: The federal Food and Drug Administration (**FDA**, United States), the International Electrotechnical Commission (**IEC**, Europe), and the Ministry of Health and Welfare (**MHW**, Japan). A GEMS MR system that is installed in a country regulated by one of these bodies has to be compliant with the limits they defined. Set the “Governing Body” value for the site as follows:
- Select the **Hardware** Tab.
 - Check/change the “Governing Body” selection so it is correct for the country the scanner is located in (**iec**, **mhw1**, **mhw2**, or **fda2**).
7. When satisfied, push the **[Configure]** button at the bottom of the tab page or wait until all tabs have been modified and go to the “**Verification**” tab. Review changes made are correct, then push the **[Configure Tab(s)]** button to update system (see Illustration 1-4).

Guided Install

File Configuration Help

Display/Add Devices Patches Reset Image/Error DB SprSnap Log file

HIS/RIS Save/Restore Option Keys SCP Service SW End Image Disk/Insite

Set Exam Number Set Disk Limit DASM VCR Setup Set Time/Date Verification

Home Network Info System Configure Hardware Configure Edit Patient Data

No Change Made **Configure tab(s)** Go to Tab

Status	Parameters	Value
	Host IP Address	3.7.25.1
	Netmask	0xffffffff
	Hostname	t1
	Domain Name	mrceg
	YP(NIS) Service	On
	IWS PC IP Address	192.168.4.1
	Magnet Monitor IP Address	1.1.1.1
	FE Laptop IP Address	192.168.4.2
	TPS Subnet	216.33.27
	Default Router IP Address	
	Power Monitor PC IP Address	
	Hospital Name	GEMS
	Suite ID	mr
	Host ID	t1
	Unique ID	1234
	Service ID	1234
	Date Format	(DD MMM YY AM/PM)
	Weight Unit	Pounds
	Monitor Type	Gray
	Palette	Gray
	Keyboard	American
	Language	English
	Field Strength	1.5T
	Package Name	ECHOSPEED SR120
	Gradient Type	8645

What's New Configure Prev Tab Next Tab

Messages

Status

Set Time Zone, Date and Time

VERIFICATION MENU – CONFIGURE TABS
ILLUSTRATION 1-4

8. Before exiting the GUI, insert an MOD into the drive with side A up. Go to the top left corner of the Install GUI and select **<File>** and then **<Save GI Configuration to MOD>**. This process **does not** create a SaveInfo disk. It just creates a copy of the information already entered in the GUI tabs for use in the next software install. Label the MOD with "Configuration Saved" include site and date.
9. Exit the install GUI by selecting **<File>** then **<Quit>**. If any error conditions still exist, you will be warned that no changes will be made before exiting.

1-2 Recovery for system with YP(NIS) Service set ON by mistake

Problem Description: User changed Network YP option to **On**. Exited GUI and initiated shutdown. Rebooted. System hung at message "setting default route". Could not proceed any further.

Solution:

1. Reset SGI computer.
2. From Startup Menu, select **[Restart]**.
3. Select **[Stop For Maintenance]**.
4. Select **Command Monitor**
5. Type: **help <Enter>** (displays list of commands)
6. Type: **single <Enter>** (for single user mode)
7. Type: **operator <Enter>** when prompted for root password.
8. Press **<Enter>** at prompt (TERM = (iris-tp) (brings up unix prompt)
9. Type: **cd /etc/config**
10. Type: **ls -al yp**
-rw-r--r-- 1 root sys <date time> yp
11. Type: **more yp** (to display contents of yp file)
*At this point the contents of the file will be either **on** or **off***
12. If contents = **on**, then type: **jot yp <Enter>** (to edit the yp file)
13. Change the contents to **off**, then Save and exit editor.
14. Reboot the system.

2- MONITOR GAMMA VALUE RESET

Do the following steps to update the monitor configurations. This insures compatibility between monitor and camera.

1. In the C-Shell, type **su <Enter>** to log in as root; enter the root password (default is **operator**).
2. **B/W Monitor:**
 - a. Type **jot /usr/g/grayPix/grayPix.cfg <Enter>** to edit the grayPix.cfg file.
 - b. Modify the numerical value for "gamma" to **0.9**

Color CRT Monitor:

- a. Type **jot /usr/g/colorPix/colorPix.cfg <Enter>** to edit the colorPix.cfg file.
- b. Modify the numerical value for "gamma" to **1.1**

Color LCD Monitor:

- a. Type **jot /usr/g/colorPix/colorPix.cfg <Enter>** to edit the colorPix.cfg file.
- b. Modify the numerical value for "gamma" to **1.05** (NEC 2010) or **1.0** (EIZO L660)

3. In the jot window, select **File -> Save**, then **File -> Exit** to save the changes. This modifies the gamma value so whenever the system reboots, the new value is not lost.
4. In the C-Shell, enter one of the following:

B/W Monitor: **gamma 0.9 <Enter>**

Color CRT Monitor: **gamma 1.1 <Enter>**

Color LCD Monitor (NEC 2010): **gamma 1.05 <Enter>**

Color LCD Monitor (EIZO L660): **gamma 1.0 <Enter>**

This is in addition to steps 2 and 3. This step immediately changes the gamma value without the need to reboot the system.


3- LOAD SERVICE SOFTWARE

Go to the install GUI window, and select the **ServiceSW/2nd Image Disk/InSite** tab and proceed to Section 3-1, Load Restricted Service or Section 3-2, Load Advanced Service. *(Restricted Service Software is only available for GE use; Advanced Service Software is only available to sites with a valid Advanced Service Package Limited License.)*

3-1 Load Restricted Service (Only available for GE use.)

To load Restricted Service Software on the system:

1. Place the appropriate Restricted Service CD-ROM in the host CD-ROM drive.
2. Select **[Load Restricted Service]**.
3. Read the proprietary message window, and respond accordingly.
4. A message "Restricted software load in progress..." will display. The load will take 5 to 10 minutes depending on host type.
5. Illustration 3-1 will display near the end of the software load. GEMS Field Engineers loading GEMS Restricted Service Software should also configure InSite at this time (see Illustration 3-1). If you do not have an InSite modem or you will configure it later, click on **[Exit]**.
6. Refer to procedure "Insite/IIP Installation" to complete the Insite configuration on the Service Methods CD-ROM.
7. Wait for the message "Restricted Service load done" in the Install GUI status bar (at bottom of window).

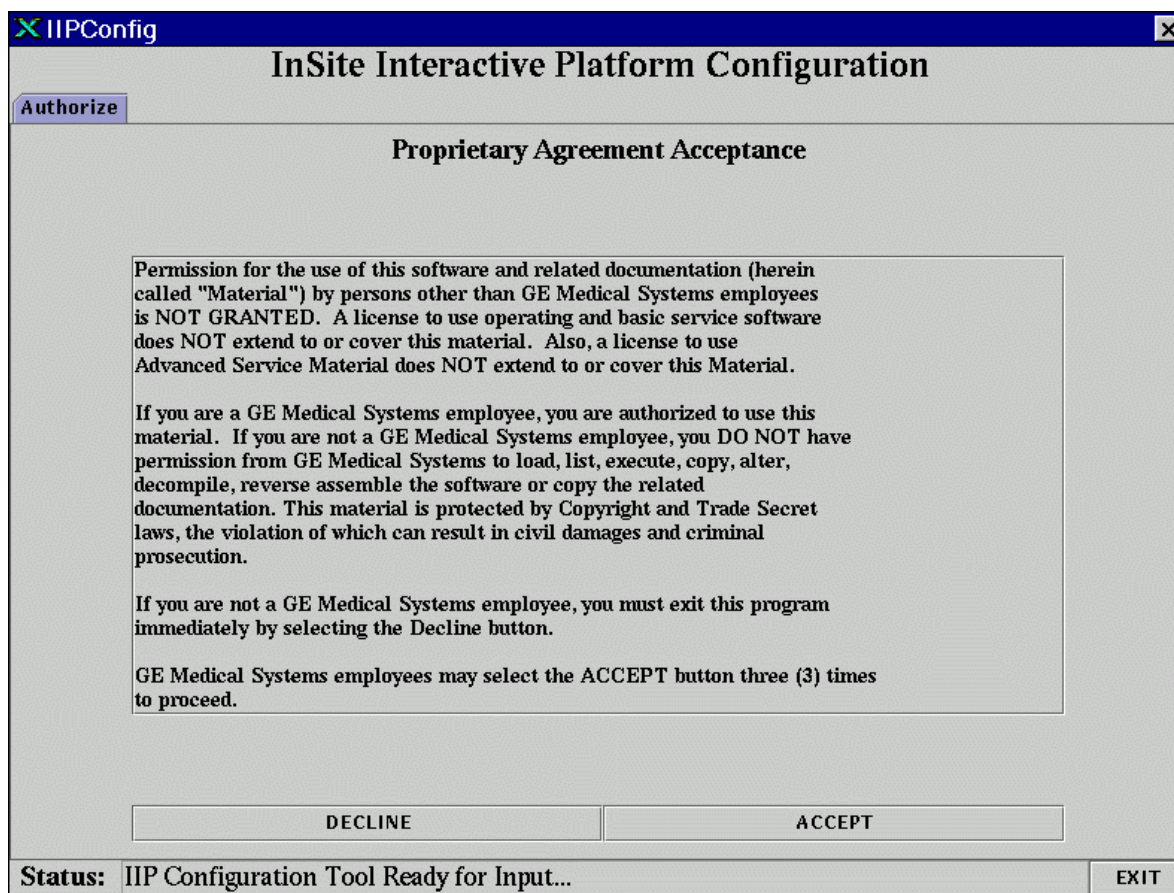

PROPRIETARY AUTHORIZATION TAB
ILLUSTRATION 3-1

3-2 Load Advanced Service (Only for sites with a valid Advanced Service Package Limited License.)

To load Advanced Service Software on the system:

1. Place the appropriate Advanced Service CD-ROM in the host CD-ROM drive.
2. Select **[Load Advanced Service]**.
3. Read the proprietary message window, and respond accordingly.
4. A message "Advance software load in progress..." will display. The load will take 5 to 10 minutes depending on host type.
5. Wait for the message "Advanced Service load done" in the Install GUI status bar (at bottom of window).

5- PERFORM SAVE INFO

After completing final site configurations, perform a SaveINFO:

1. Put the MOD that you previously used to **Save GI Configuration to MOD** with side A up (same side is OK because a separate partition will be created for the SaveINFO).
2. Select the **Save/Restore** tab in the Install GUI and Save the system info.

Note

If the MOD has not been used, you will see a message regarding the use of the MOD for “Raw Data Storage”; choose **[Overwrite]**.

3. SaveInfo is now complete. Remove the InfoMOD and write the date and time on it.

6- FINAL STEPS

1. Exit the install GUI by selecting **<File>** then **<Quit>**.
2. On the Service Desktop, click on **[System Shutdown]**.
3. After the system is all the way down, when prompted, click on **[Restart]**.
4. On the login screen, login as **Signa** and bring the system up for the configuration changes to take effect.
5. Complete any remaining system calibrations per the installation/upgrade flowchart.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	Mar. 10, 2000	J. Wolak	Initial version document.
1	Sept. 21, 2000	M. Keber	Updated to include Release 8.4 and per SPR MRlge60632 to add "Governing Body" check for new installations.
2	Dec. 14, 2000	D. Hofstetter	Moved "Governing Body" check to Section 1. Updated all graphics to show correct tab locations for 8.4 W/ACGD, ASP1 and ASP2 Releases.
3	Aug. 13, 2001	J. Wolak	Added Release 9.0 to note before Sec. 1. Added gamma settings for NEC 2010 & EIZO L660 monitors in Sec. 2