

TABLE OF CONTENTS

TABLE OF CONTENTS.....	1
1- TOOLS REQUIRED.....	2
2- REMOVAL OF FAILED RF/PEN I/F PANEL.....	2
3- INSTALLATION OF NEW I/F PANEL	3
4- FUNCTIONAL CHECKS REQUIRED.....	4
REVISION HISTORY	5

Description - This section applies only to Signa 1.5T or 1.0T systems installed with the Erbtec RF/Pen cabinet.

1- TOOLS REQUIRED


#2 nonferrous Phillips screwdriver

WARNING!


POSSIBLE PERSONAL INJURY! REMOVE POWER AND ALL EXTERNAL CABLES FROM THE RF/PEN CABINET BEFORE YOU REMOVE ANY COMPONENTS FROM IT. VERIFY THAT ALL BREAKERS IN THE REAR OF CABINET ARE OFF, AND DISCONNECT MAIN POWER FEED CORD BEFORE OPENING CHASSIS FOR SERVICE. SERIOUS INJURY OR DEATH BY ELECTROCUTION MAY OTHERWISE OCCUR. LOCK OUT AND TAG OUT THE RF/PEN CABINET.

2- REMOVAL OF FAILED RF/PEN I/F PANEL

1. Turn off RF/Pen cabinet power circuit breakers on rear of RF/Pen cabinet. Lock out and tag out the RF/Pen cabinet. Illustration L3129a shows the location of the RF/Pen interface (I/F) panel on the back of the RF/Pen cabinet. (Refer to *Procedure For Safety: Section 6*)

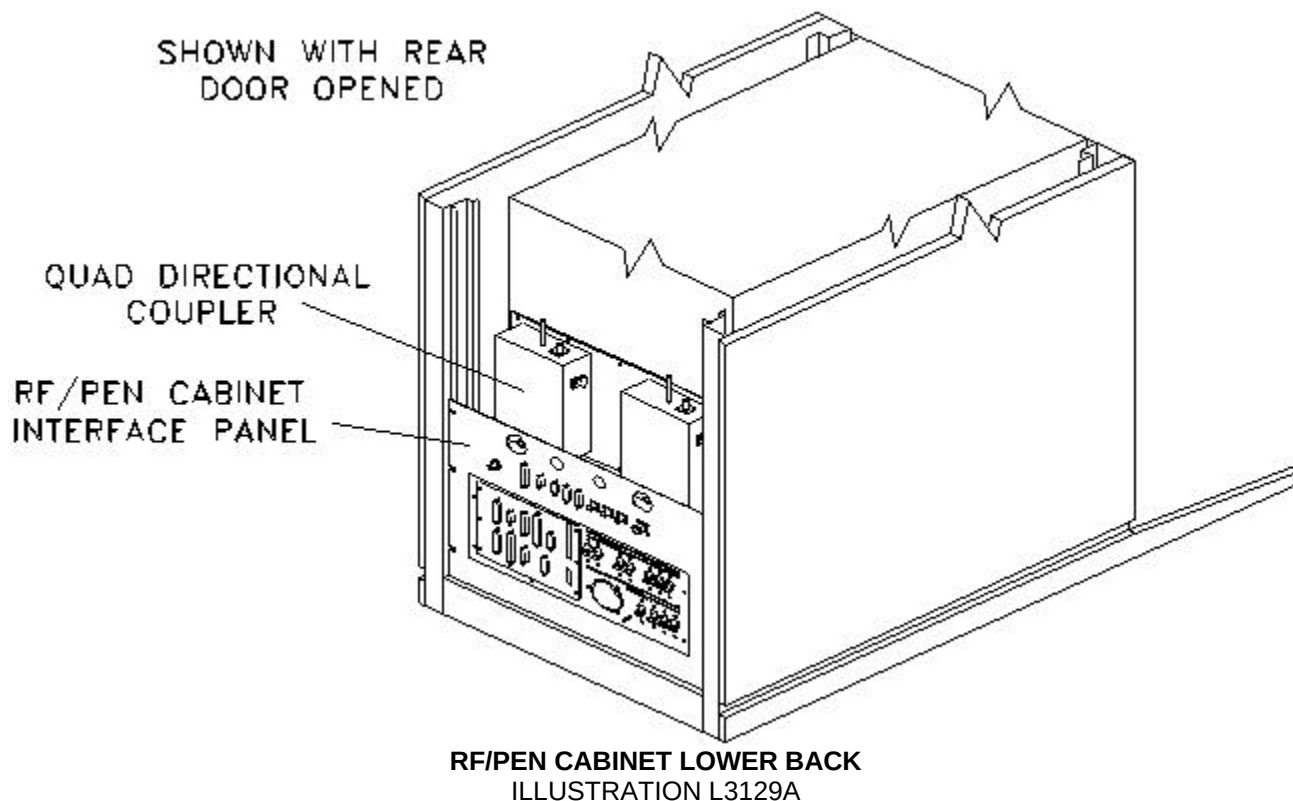

2. Verify that power to the RF/Pen Cabinet is off by checking the fans and LEDs on both the RF amplifier and the RF System Controller (RFSC).
3. Disconnect the cables on the I/F panel and on the rear of the magnet enclosure power supply (MEPS) module.
4. Disconnect the optical cable and RF power cables that run through the I/F panel.
5. Remove the six 10-32 X 1/2 screws and six #10 flat washers securing the I/F panel to the RF/Pen cabinet. See Illustration L3141a.

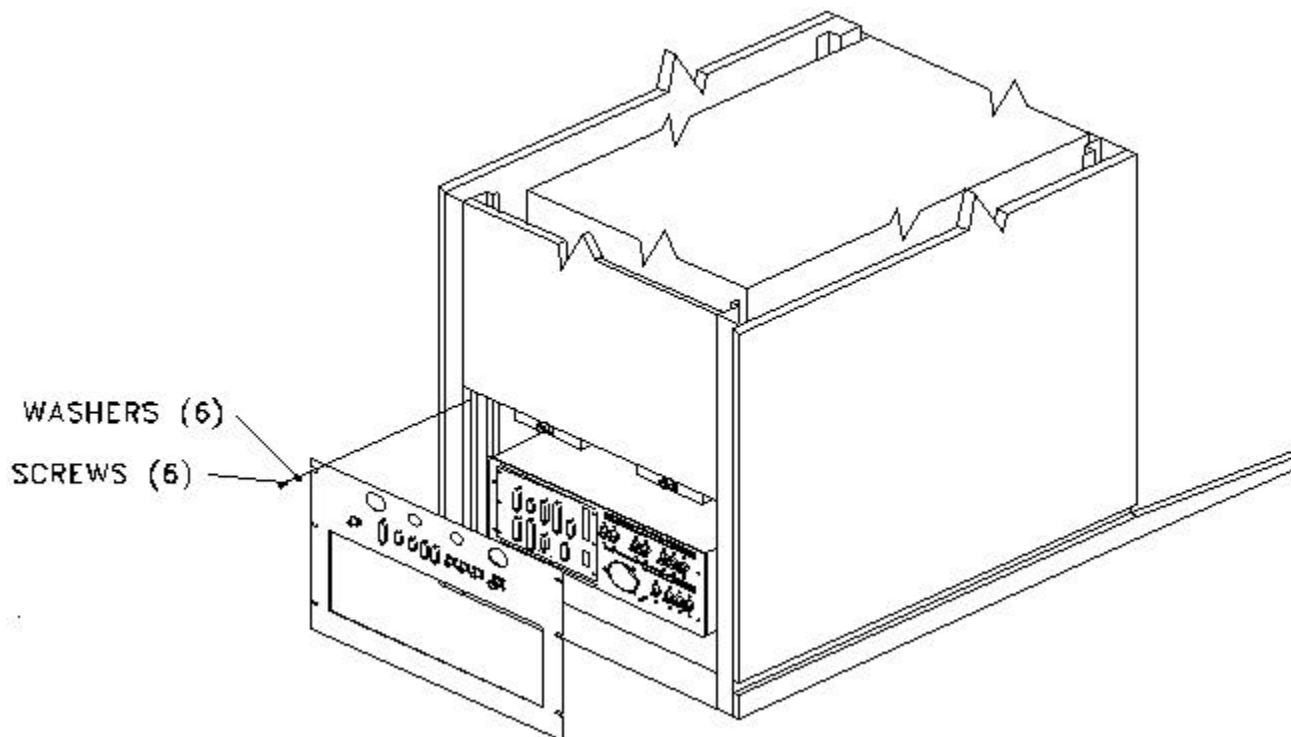

RF/PEN CABINET INTERFACE PANEL REMOVAL
ILLUSTRATION L3141A

6. Remove the I/F panel.

3- INSTALLATION OF NEW I/F PANEL

1. Position the I/F panel on rear of RF/Pen cabinet.
2. Secure with six 10-32 X 1/2 screws and six #10 flat washers. See Illustration L3141A in Section 2.
3. Connect the cables to the I/F panel and the MEPS module, matching the labels on the cables to the appropriate connection.
4. Connect the optical cable and RF cables through the I/F panel.

5. Remove the lock out/tag out, and restore power to the RF/Pen cabinet. There are no adjustments or calibrations to be made.

4- FUNCTIONAL CHECKS REQUIRED

1. Perform one head or body scan.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	July 15, 1998	J. Saperstein	Initial conversion from Toolbook to Word.
1	May 21, 1999	SM Atladottir	Updated Procedure References for New GUI