

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1- 1.5T BODY T/R SW & HYBRID ASM., IDD, AND PREAMP	2
1-1 Removing Old 1.5T Body T/R Hybrid Hardware.....	2
1-2 Installing New 1.5T Body T/R Hybrid Hardware.	3
1-3 Required Calibrations/Checks.....	5
2- COMBINED 1.5T OR 1.0T BODY T/R HYBRID HARDWARE.	6
2-1 Removing Body T/R Sw, & Hybrid Asm. and Preamp	6
2-2 Installing Body T/R Sw, & Hybrid Asm. and Preamp	7
2-3 Required Checks.....	8
REVISION HISTORY.....	9

This procedure applies to 1.5T and 1.0T systems.

Description - Procedure for replacing the Body TR Switch and Preamplifier assemblies, which are located in the rear pedestal.

1- 1.5T BODY T/R SW & HYBRID ASM., IDD, AND PREAMP

Note

This section applies to Signa Advantage 1.5T system only. In 1993, the Inductive Drive Disable Switch Boxes were replaced with an Inductive Drive Disable Board which is located in the Body T/R Switch and Hybrid Assembly (46-321312G1). This new assembly is the FRU for the old Inductive Drive Disable Switch Boxes, old Body T/R Switch and Hybrid Assembly, or the new Body T/R Switch, IDD, and Hybrid Assembly.

1-1 Removing Old 1.5T Body T/R Hybrid Hardware



Possible equipment damage. Verify that power is removed from RF cabinet before removing or disconnecting exam room RF components. Verify that the breaker on rear of cabinet is off before starting work. The RF amplifier or TR driver could be damaged by operating into improperly terminated loads.

1. Turn off RF cabinet power circuit breaker, CB1, on rear of RF cabinet. Lock out and/or tag out cabinet. (Refer to CD-ROM *Dir. 2187583-3 [or -2], MR Release Signa 5x/8x Service Methods, Renewal Parts and Service Tools*, Procedure for Safety, Section 6, OSHA LOCKOUT/TAGOUT REQUIREMENTS.)
2. Verify that power to the RF/Pen II Cabinet is off by checking the fans and LEDs on both the RF amplifier and the System Support Module (SSM).
3. Remove the side covers on the rear pedestal.
4. Disconnect the cables from the Body T/R Switch & Hybrid Assembly, Preamp and Drive Disable Assembly. See Illustration L2921A .

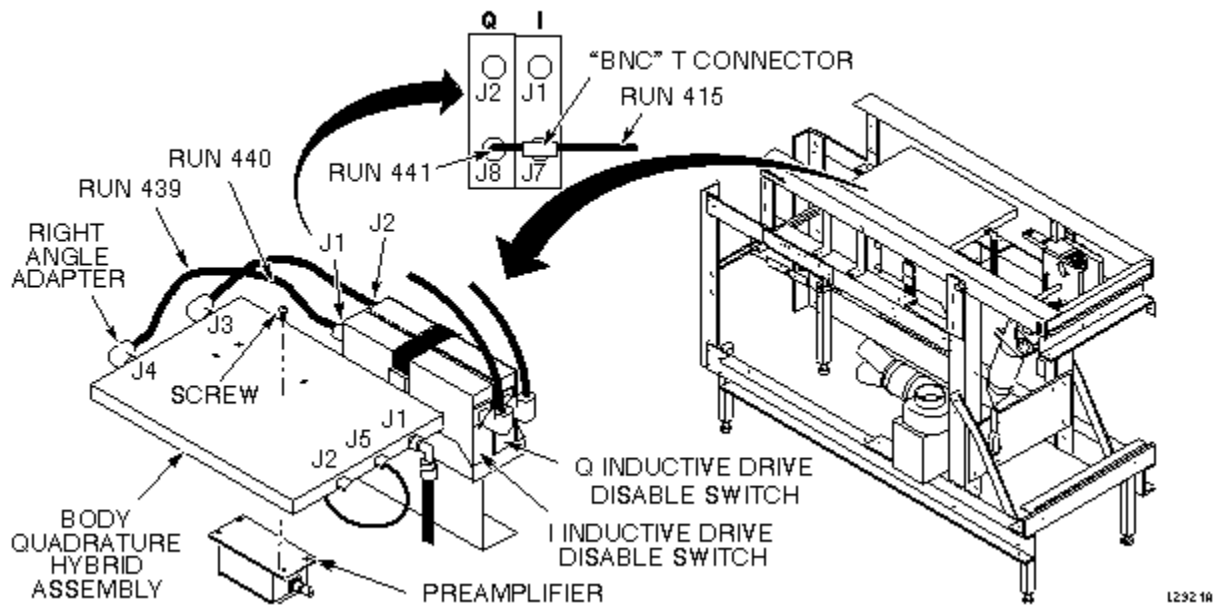


ILLUSTRATION L2921A
1.5T BODY T/R SWITCH & HYBRID AND PREAMP

5. Remove the Body T/R Switch & Hybrid Assembly, Preamp and Drive Disable Assembly from the rear pedestal.
6. Remove the Preamp and Drive Disable Assembly from the Body T/R Switch & Hybrid Assembly.

1-2 Installing New 1.5T Body T/R Hybrid Hardware.

Note

This installation procedure is for the new Body T/R Switch, IDD, and Hybrid assembly which is the FRU for the old Inductive Drive Disable Switch Boxes or old Body T/R Switch and Hybrid Assembly.

1. Remove the preamplifier and coaxial connector from the original Body Hybrid (46-264170G1), and install it on the new combined Body T/R Switch & Hybrid - IDD assembly (46-321312G1), exactly as they were installed on the original. See Illustration L2922A .

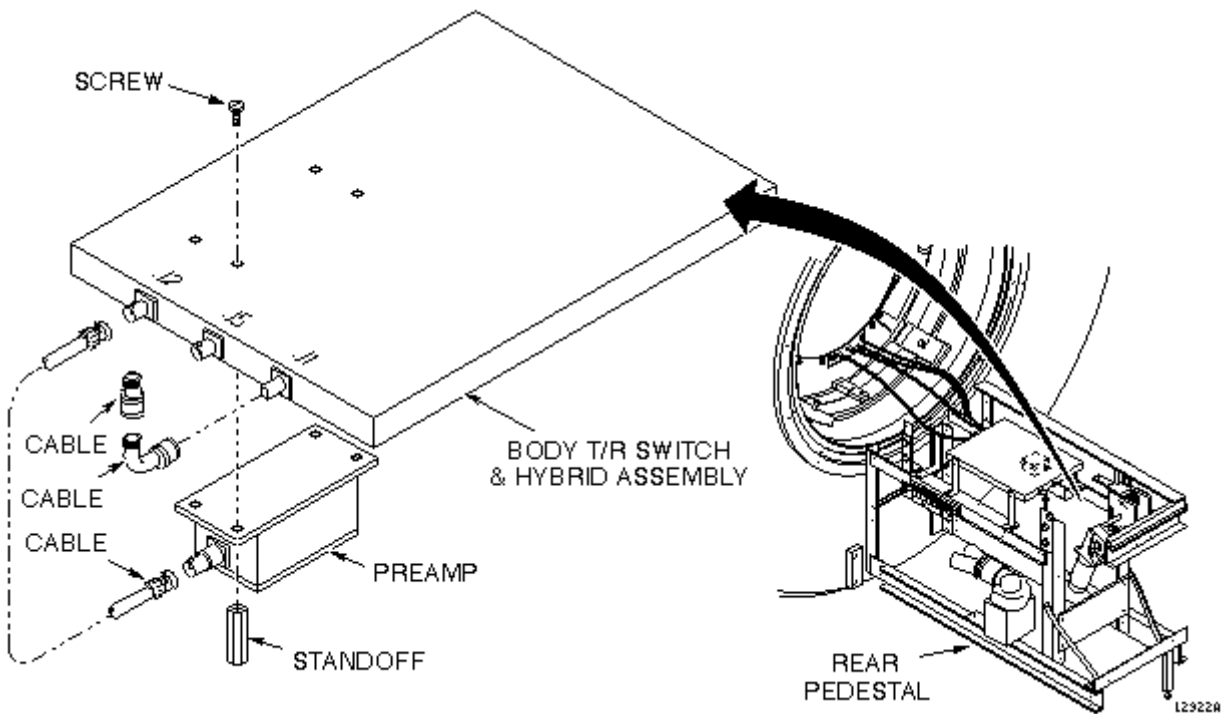


ILLUSTRATION L2922A
1.5T & 1.0T BODY T/R SWITCH & HYBRID AND PREAMP

2. Install the replacement Body Hybrid assembly on the same brackets that the old assembly was removed from.
3. The Inductive Drive Disable boxes and the two 1/4 wave length cables (46-243775G682 and 46-243775G683) will not be reused. The Inductive Drive Disable boxes should be returned with the Body T/R Switch & Hybrid assembly to close out the exchange. See Illustration L2921A .

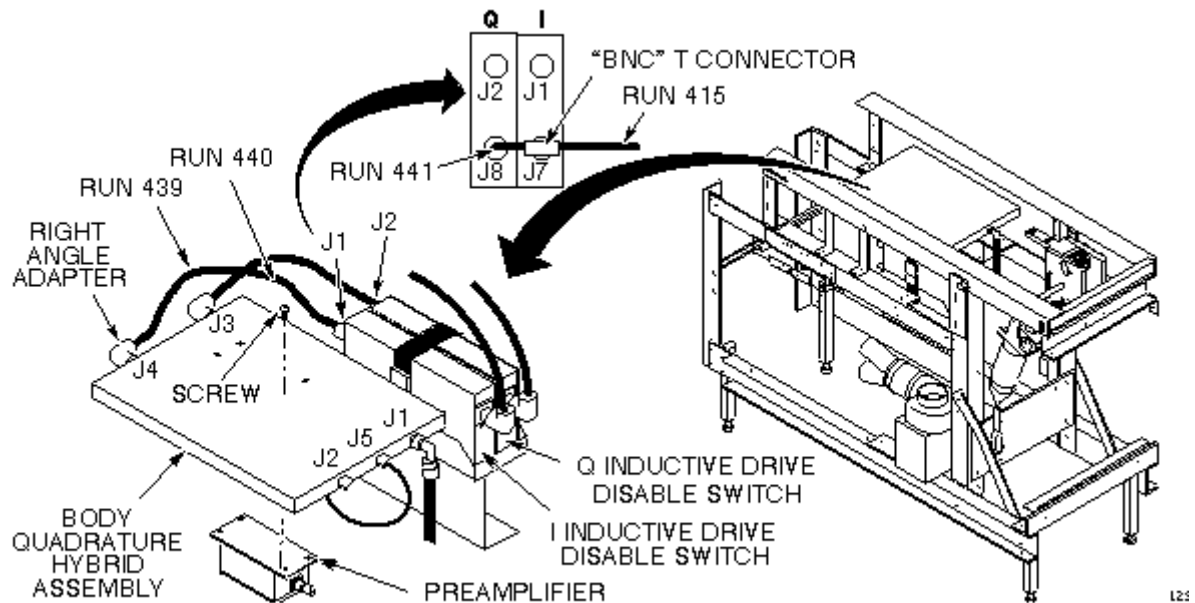


ILLUSTRATION L2921A
1.5T BODY T/R SWITCH & HYBRID AND PREAMP

4. Remove the coaxial T adapter, and short coaxial cable (run 441) from the Inductive Drive Disable bias line (run 415), that were formerly connected to MG3A34J7 and MG3A33J8. Connect this line to J6 of the new combined Hybrid Inductive Dive Disable (46-321312G1). Change the cable label to read MG3-A3-J6.
5. Connect the 1/2 wavelength cables (46-287688G1 & 46-287688G2), that were disconnected from the Inductive Drive Disable boxes to the new hybrid assembly on connectors J3 and J4. The opposite end of these cables should still be connected to the body coil at J1 and J2. If the system has a normal GE magnetic field polarity, then the cable from Body Coil J1 connects to MG3-A3-J4 (Hybrid J4), and Body Coil J2 connects to MG3-A3-J3 (Hybrid J3). If the systems Magnet is reverse ramped, the Body Hybrid connections previously stated should be reversed. Relabel these cables accordingly.
6. Connect the RF Body Heliax to the Hybrid connector J1. Connect the Preamplifier out cable (run 262), to the body preamplifier output connector. The connector J5 is the service port and has nothing connected to it.
7. Install the side covers on the rear pedestal.
8. Remove lock out and/or tag out from the RF cabinet. Turn on RF cabinet power circuit breaker, CB1, on rear of RF cabinet.

1-3 Required Calibrations/Checks

1. Calibration is required when replacing the old 1.5T Body Hybrid hardware with the new version hardware. Refer to procedure for DD/TR Driver Board Adjustments.

2. Check body SNR. Refer to procedure for Signal To Noise Check, Section 2- Body Scans.

(An alternative proprietary procedure is available for GE use, or to customers with an Advanced Service Package Limited License. Refer to procedure for TLT and perform a Body SNR scan.)

2- COMBINED 1.5T OR 1.0T BODY T/R HYBRID HARDWARE.

This procedure applies to the new combined Body TR Switch & Hybrid Assembly with Preamp. See Illustration L2922A .

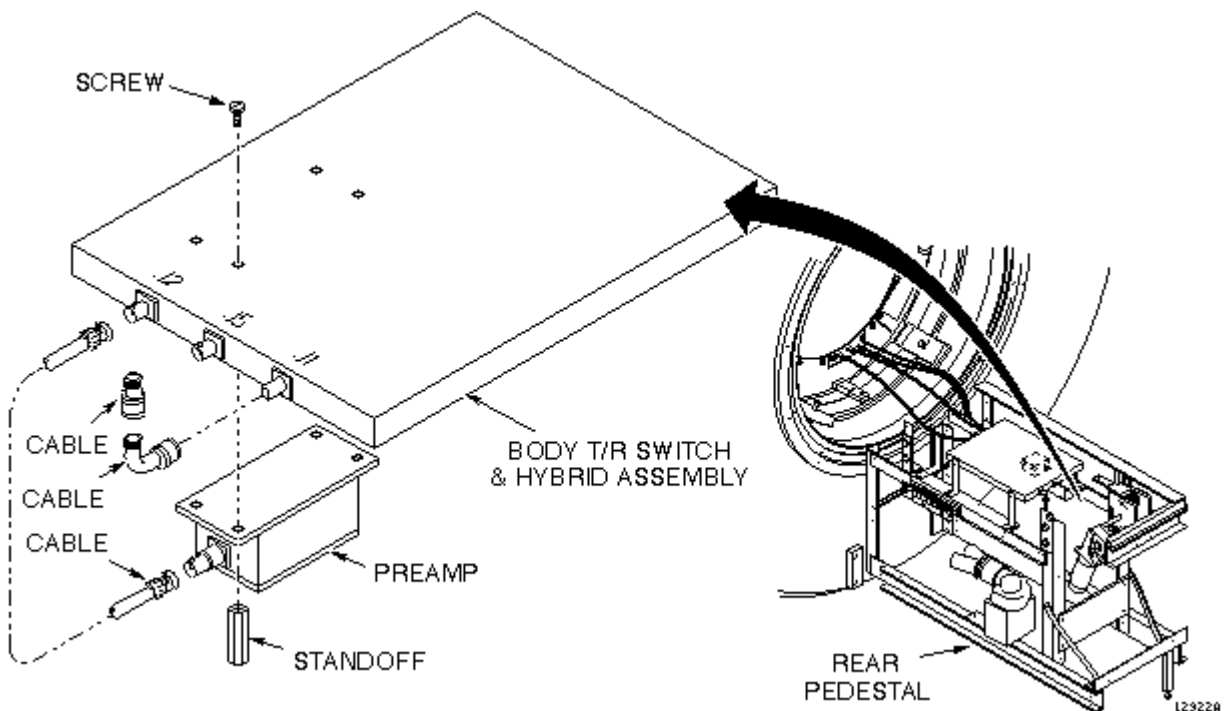


ILLUSTRATION L2922A
1.5T & 1.0T BODY T/R SWITCH & HYBRID AND PREAMP

2-1 Removing Body T/R Sw, & Hybrid Asm. and Preamp

CAUTION

Possible equipment damage. Verify that power is removed from RF cabinet before removing or disconnecting exam room RF components. Verify that the breaker on rear of cabinet is off before starting work. The RF amplifier or TR driver could be damaged by operating into improperly terminated loads.

1. Turn off RF Cabinet power circuit breaker, CB1, on rear of RF Cabinet. Lock out and/or tag out the RF Cabinet. (Refer to CD-ROM *Dir. 2187583-3 [or -2], MR Release Signa 5x/8x Service Methods, Renewal Parts and Service Tools*, Procedure for Safety, Section 6, OSHA LOCKOUT/TAGOUT REQUIREMENTS.)
2. Verify that power to the RF/Pen II Cabinet is off by checking the fans and LEDs on both the RF amplifier and the System Support Module (SSM).
3. Remove the side covers on the rear pedestal.
4. Disconnect the cables from the Body T/R Switch & Hybrid Assembly and the Preamp. See Illustration L2922A .

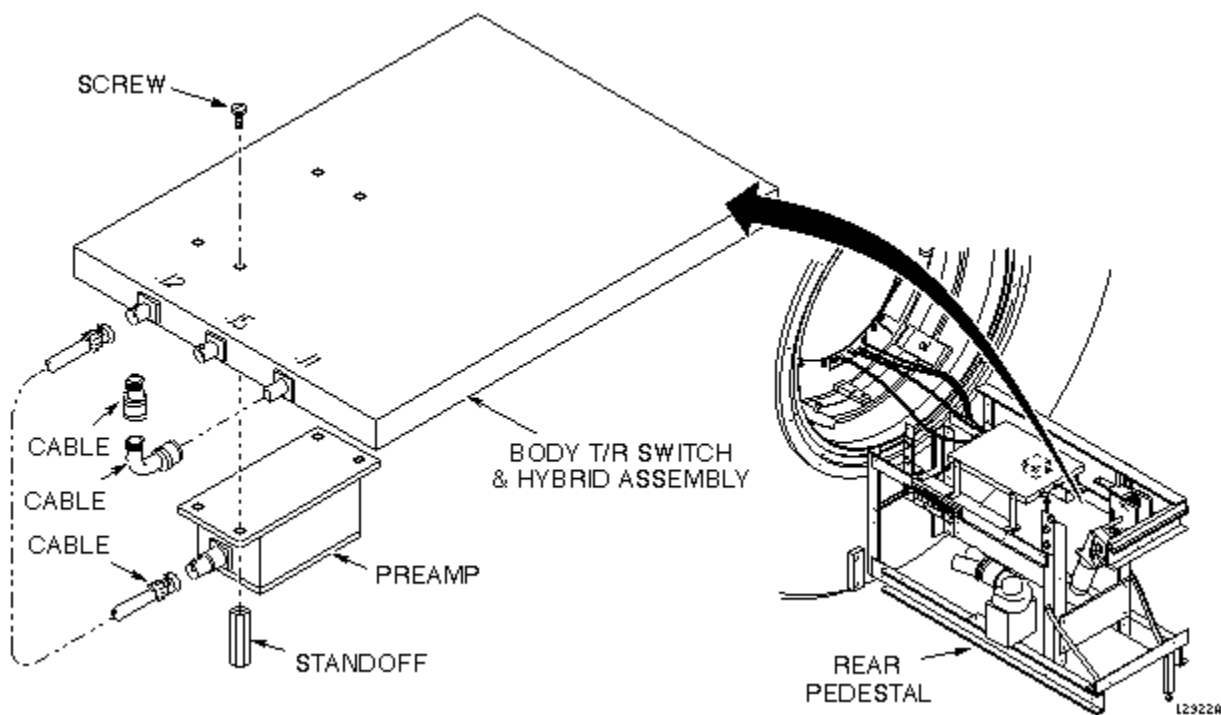


ILLUSTRATION L2922A
1.5T & 1.0T BODY T/R SWITCH & HYBRID AND PREAMP

5. Remove the Body T/R Switch & Hybrid Assembly and the Preamp from the rear pedestal.
6. Remove the Preamp from the Body T/R Switch & Hybrid Assembly.

2-2 Installing Body T/R Sw, & Hybrid Asm. and Preamp

1. Refer to Section 2-1 for removal of the Body T/R & Hybrid Assembly and preamp.
2. Install the preamp to the Body T/R Switch & Hybrid Assembly. See Illustration L2922A.
3. Install the Body T/R Switch & Hybrid Assembly and the preamp to the rear pedestal.

4. Connect the cables to the Body T/R Switch & Hybrid Assembly and the Preamp. See Illustration L2922A.
5. Install the side covers on the rear pedestal.
6. Remove lock out and/or tag out from the RF Cabinet. Turn on RF Cabinet power circuit breaker, CB1, on rear of cabinet.

2-3 Required Checks

1. Refer to procedure for DD/TR Driver And Pin Diodes Checks.
2. Check Body SNR. Refer to procedure for Signal To Noise Check.

(An alternative proprietary procedure is available for GE use, or to customers with an Advanced Service Package Limited License. Refer to procedure for TLT and perform a Body SNR scan.)

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	Aug 19, 1998	R. Hawthorne	Initial conversion to Word