

TABLE OF CONTENTS

TABLE OF CONTENTS 1

1- CERD POWER SUPPLY 2

 1-1 Removal of Failed Power Supply..... 2

 1-2 Installation of Replacement Power Supply..... 3

 1-3 Functional Checks Required..... 4


REVISION HISTORY..... 5

Description - This material describes replacement of the CERD (combined exciter/receiver DAB) Power Supply, located in the ISE/TPS Chassis.

1- CERD POWER SUPPLY


1-1 Removal of Failed Power Supply

1. Remove the System Cabinet front cover.
2. Turn off the main circuit breaker, CB1 (see Illustration L4601A).

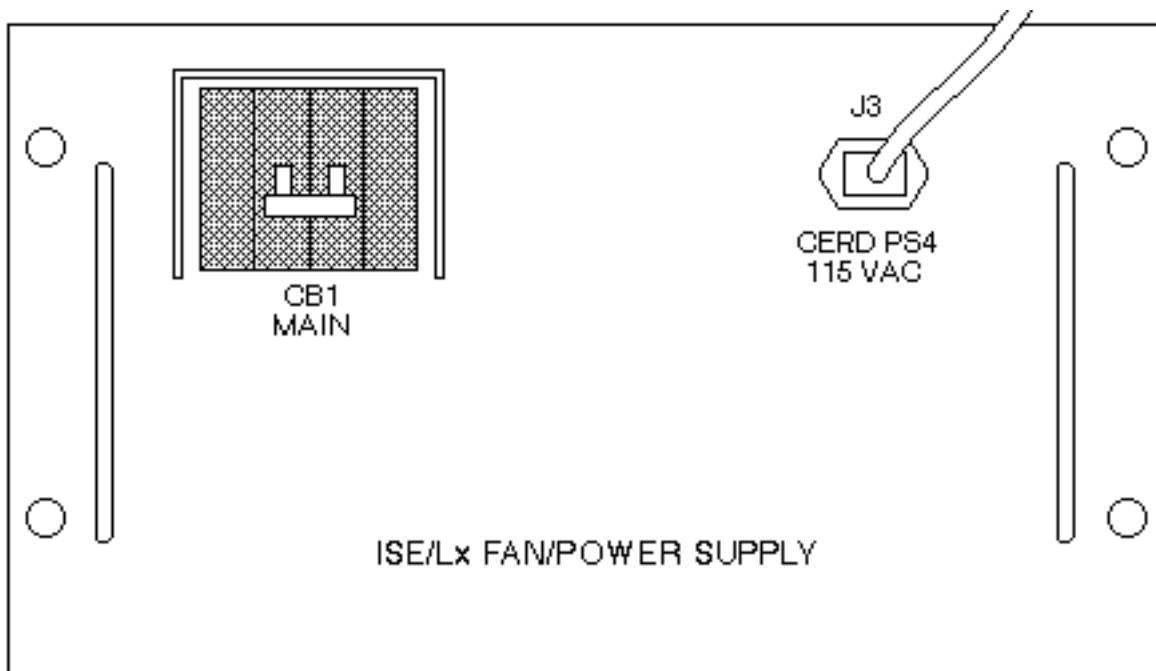

ISE/TPS FAN/POWER ASSEMBLY - FRONT VIEW
ILLUSTRATION L4601A

3. Remove cables connected to J1, J2, and J3 on the CERD Power Supply Assembly.
4. Loosen the six front panel screws. Remove the CERD Power Supply Assembly from the ISE/TPS Chassis by pulling outward with the handle that is provided on the CERD Power Supply (see Illustration L4729A).

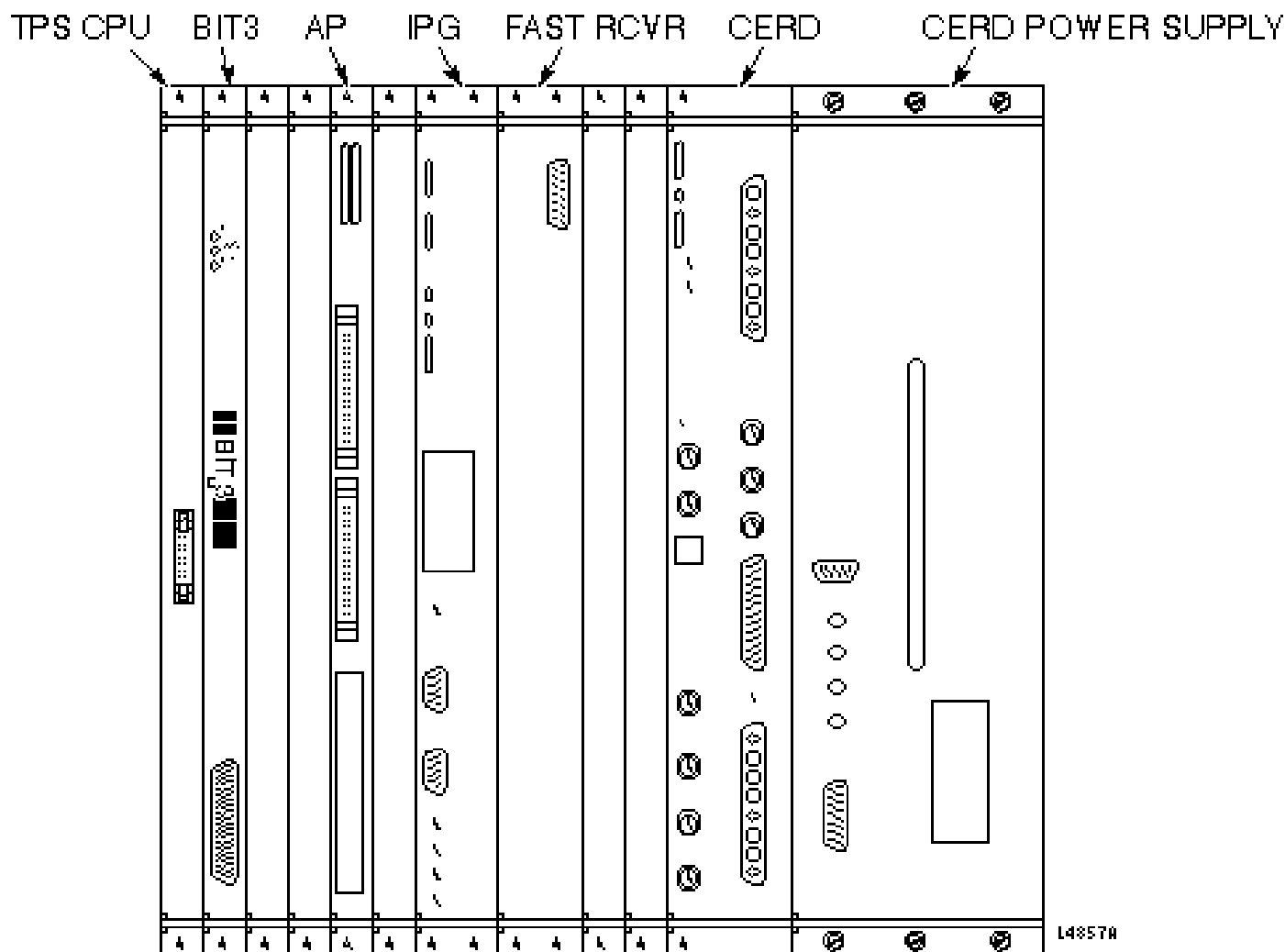

ISE/TPS CHASSIS - FRONT VIEW
ILLUSTRATION L4729A

1-2 Installation of Replacement Power Supply.

1. Install the CERD Power Supply Assembly in the ISE/TPS Chassis and tighten the six front panel screws (see Illustration L4729A in Section 1-1).
2. Reconnect cables back to J1, J2, and J3 on the CERD Power Supply Assy (see Illustration L4729A in Section 1-1).
3. Turn on the main circuit breaker, CB 1 (see Illustration L4601A in Section 1-1).
4. Refer to Procedure CERD Power Supply Adjustments for adjustments and checks of the voltages on the power supply.
5. Install the System Cabinet front cover.
6. If Signa Horizon (8.X) software is running, **[Reset TPS]**.

1-3 Functional Checks Required

1. Perform one head or body scan.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	June 1, 1998	J. Saperstein	Initial Conversion from Toolbook to Word.