

TABLE OF CONTENTS

TABLE OF CONTENTS ..... 1

1- INLINK HARDWARE OVERVIEW ..... 2

2- INLINK PREINSTALLATION ..... 2

3- INLINK HARDWARE INSTALLATION ..... 2

    3-1 Contents of InLink Kit ..... 2

    3-2 Required Tools, Supplies, And Test Equipment..... 3

    3-3 Modem And Modem Power Supply Installation..... 3

    3-4 Cable Installation..... 4

    3-5 Final Steps ..... 5

REVISION HISTORY..... 7

## 1- INLINK HARDWARE OVERVIEW

This procedure tells you how to install the InLink Kit option on an MR Signa Advantage 1.5T, 1.0T, or 0.5T (Release 5.x) System. The major tasks consist of hardware installation of the modem and cabling.

The InLink kit provides remote access capability to the Signa Advantage MR System. The equipment installed at the Customer's site consists of a 19.2 kilobaud modem connected to the system computer through an RS-232 serial cable. The network is completed by connecting the modem to a voice grade DID, Direct Inward Dial, phone line.

## 2- INLINK PREINSTALLATION

Perform the following items prior to InLink hardware installation:

- 1. General Electric Medical Systems Personnel only:** General Electric Company must have a valid InSite Package Agreement or a MSA agreement (Revision 4/92 or later) signed by the Customer and General Electric Company. Your local Area Manager should have a Customer Introduction Package for InSite.
2. Assist the Customer with locating the DID, Direct Inward Dial, voice grade line receptacle. It should be of the RJ-11 type and preferably located at the rear of the System Cabinet (MR2) location. It is the Customer's responsibility to have this line properly installed.
3. InLink kit 46-301708G7 (U. S.) ordered and received or 46-301708G8 (International) and locally procured modem and modem power supply ordered and received.

## 3- INLINK HARDWARE INSTALLATION

### 3-1 Contents of InLink Kit

InLink Kit 46-301708G7 consists of:

| QTY. | PART NUMBER   | DESCRIPTION                                                      |
|------|---------------|------------------------------------------------------------------|
| 1    | 46-278761P2   | High Speed, Dial-Up Modem with Modem Power Supply                |
| 1    | 46-278763P1   | 25' Phone Line Cord with RJ-11 Plugs                             |
| 1    | 46-251182P322 | Rating Plate                                                     |
| 1    | 46-287683P2   | "Property of GE" Label                                           |
| 1    | 46-297005G3   | 6' Cable (Modem to MR2A11J30)                                    |
| 1    | 46-271601G20  | 15' Cable (MR2A11J30 to Sun CPU Board Serial Port B)             |
| 1    | 46-136246P16  | 2' Adhesive Backed 1" Wide Loop Strip                            |
| 1    | 46-136246P18  | 2' Adhesive Backed 1" Wide Hook Strip                            |
| 2    | 46-265067P1   | Panel Interface Clip                                             |
| 1    | 46-015412     | Signa Advantage 1.5T, 1.0T, & 0.5T InLink Installation Direction |
| 1    | 46-320095P1   | InSite Cabinet Magnet - Computer (Main)                          |
| 1    | 46-320095P4   | InSite Cabinet Magnet - PDU                                      |
| 1    | 46-320095P5   | InSite Cabinet Magnet - Systems                                  |

|   |             |                                     |
|---|-------------|-------------------------------------|
| 1 | 46-320095P6 | InSite Cabinet Magnet - Penetration |
| 1 | 46-320095P7 | InSite Cabinet Magnet - Gradient    |
| 1 | 46-320095P8 | InSite Cabinet Magnet - RF          |

**Note**

Modem and modem power supply must be locally procured for InLink Kit 46-301708G8.


**3-2 Required Tools, Supplies, And Test Equipment**

Standard hand tools are used for the hardware installation. No special tools or test equipment are required.

**3-3 Modem And Modem Power Supply Installation****Note**

This procedure was written for InLink Kit 46-301708G7. The locally procured modem for InLink Kit 46-301708G8 may be different.

1. Remove System Cabinet (MR2) front cover.
2. Remove and retain screws securing bottom panel on System Cabinet.
3. Arrange existing cables in bottom of System Cabinet so there's room for the modem and modem power supply.
4. Install Modem and Modem Power Supply in System Cabinet per Illustration L1970A.

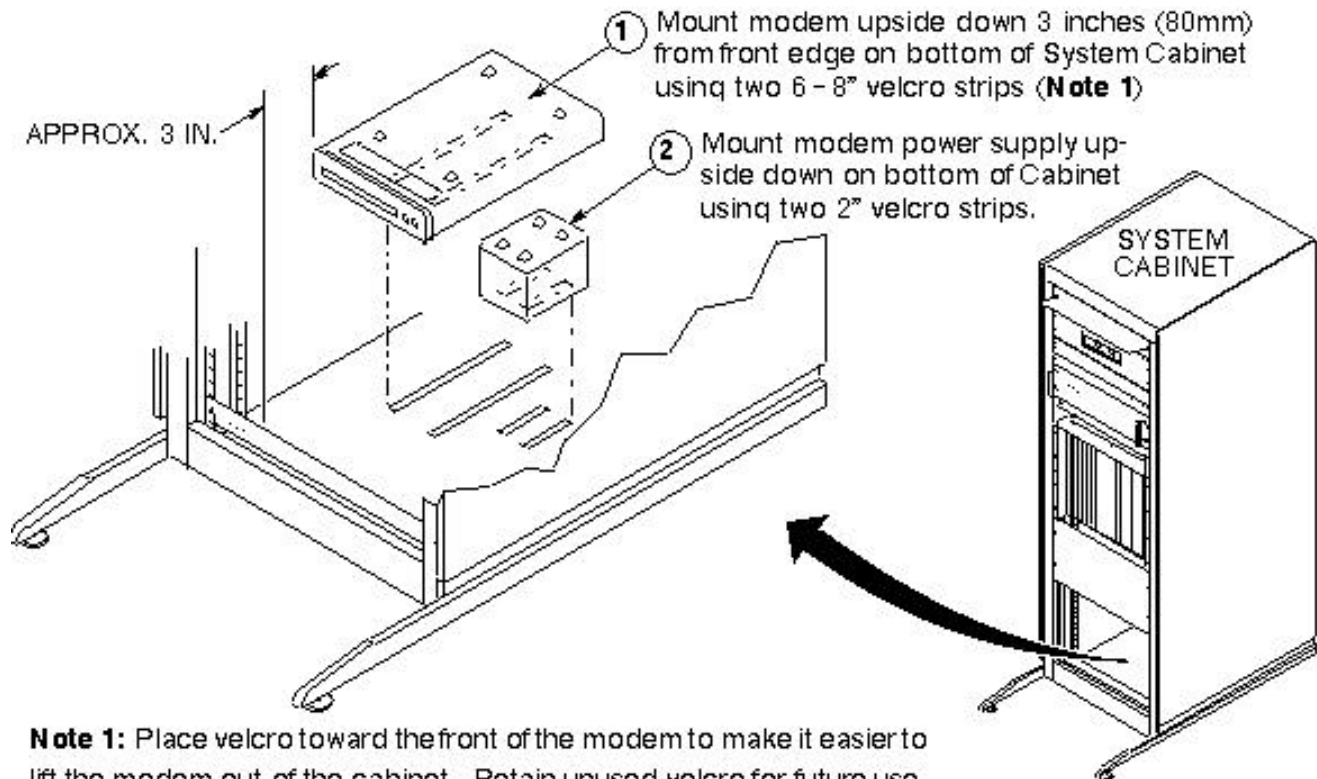

**Note 1:** Place velcro toward the front of the modem to make it easier to lift the modem out of the cabinet. Retain unused velcro for future use.

ILLUSTRATION L1970A  
**MODEM INSTALLATION IN SYSTEM CABINET**

### 3-4 Cable Installation

#### **Note**

This instruction was written for InLink Kit 46-301708G7. The locally procured modem for InLink Kit 46-301708G8 may be different.

Install InLink hardware cabling. Modem and modem power supply can be lifted as necessary to facilitate cable installation. See Illustration L1975A.

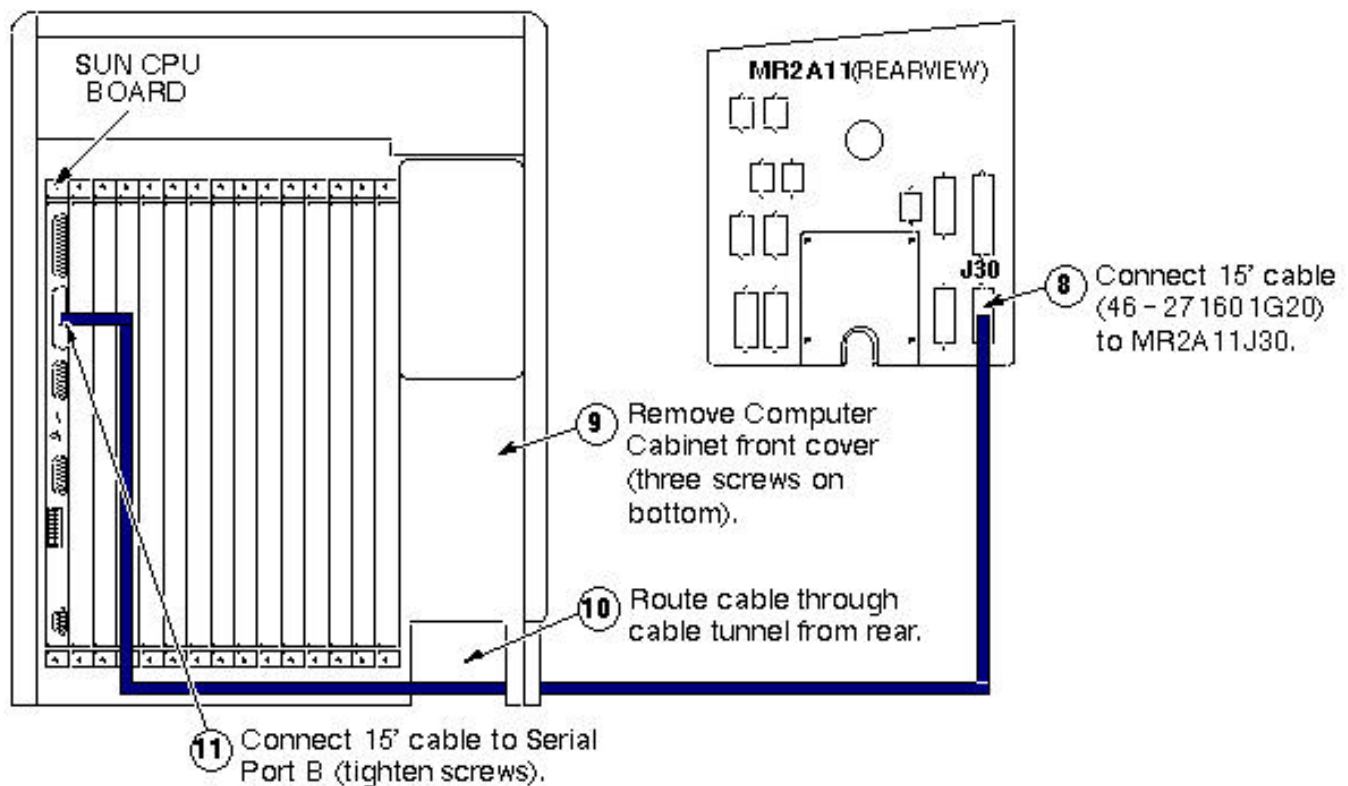
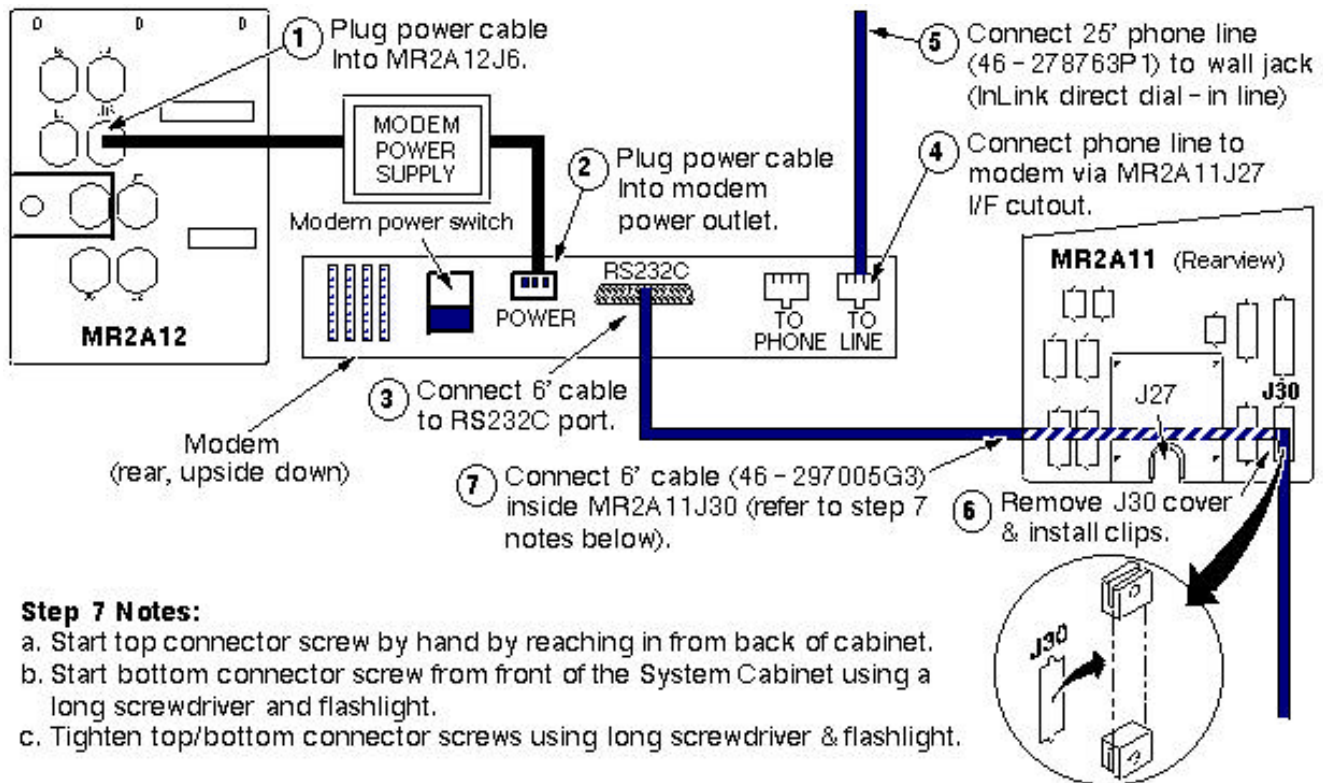


ILLUSTRATION L1975A  
INSITE HARDWARE INTERCONNECTIONS

### 3-5 Final Steps

1. Install Computer Cabinet front cover.

2. Ensure modem and modem power supply are securely mounted in System Cabinet.
3. Install System Cabinet bottom panel and front cover.
4. Install rating plate along top edge on rear of System Cabinet.
5. Install ID Magnets on front of appropriate cabinets. Labels can be placed on either right or left upper corner, depending on customer preference.

**Note**

Use adhesive-back tape or left over velcro to install the Computer magnet onto its non-magnetic cabinet.

6. **For GE Medical Systems personnel only:** Load MR InSite software and perform InSite notification/checkout (refer to Procedure for Load/Delete InSite Software).

## REVISION HISTORY

| REV | DATE           | AUTHOR     | PRIMARY REASONS FOR CHANGE                |
|-----|----------------|------------|-------------------------------------------|
| 0   | August 5, 1998 | M. Whitlow | Initial conversion from Toolbook to Word. |
|     |                |            |                                           |
|     |                |            |                                           |
|     |                |            |                                           |
|     |                |            |                                           |