

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1- INTRODUCTION	2
1-1 Procedure.....	2
REVISION HISTORY	6

1- INTRODUCTION

Image Snap is a tool used to aid in troubleshooting image quality issues. An image can be captured and stored in a data file for automatic transfer and review to the Automated Support Center (ASC). The following will be captured into a file:


- Image data files
- Raw data files (if they exist)
- System operating information
- Comments entered by operator

Note

Insite is required for auto transfer to occur. However, the files will be saved if Insite is not available and/or there is a failure during auto transfer.

1-1 Procedure

Image Snap is accessed by selection from the Service Desktop [**Service Desktop**] then [**ImageSnap**]. (see Illustration 1-1)

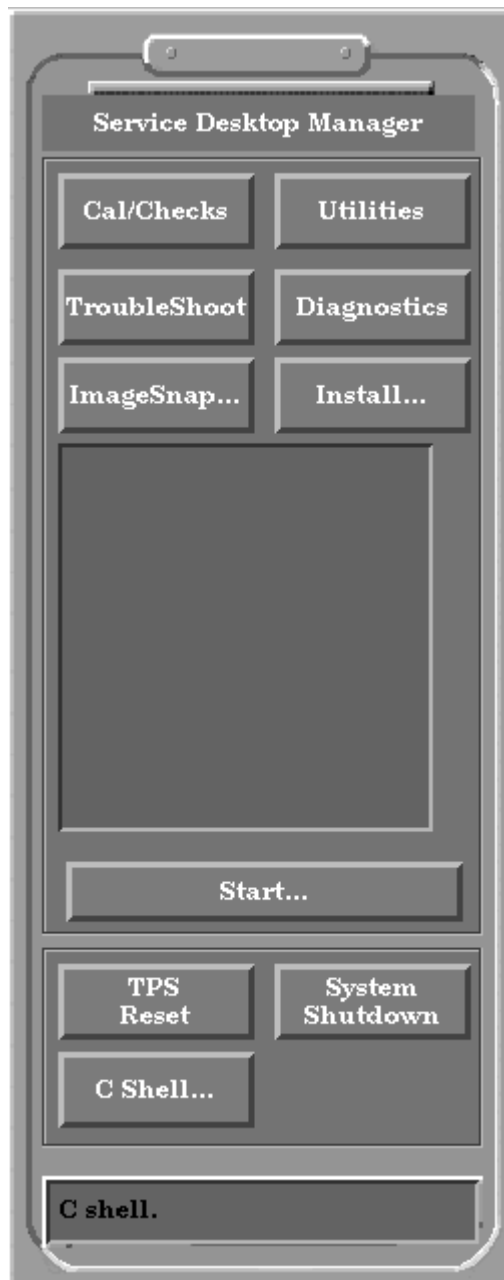

SELECTING IMAGE SNAP
ILLUSTRATION 1-1

After selecting **[ImageSnap]** a screen opens with available exams and series on the system disk. The most current exam and series is automatically selected. If another exam or series is desired, click on that exam or series. If Raw Data exists, the Yes box will be highlighted. Place cursor in the Comment field to enter comments; files will not be sent without an entry in the Comment field. Up to 3 image files can be selected for the exam/series selected. A dialog box will be displayed if more than 3 image files are selected. When the exam or series selection has been saved, select **[File ...]** to verify that the file has been saved. (see Illustration 1-2)

Note

The **[Save]** button will not be selectable until an image is highlighted.

The screenshot shows a window titled "Image Snap Shot" with three columns: "Exam:", "Series:", and "Image:". The "Exam:" column contains a list of exams with IDs and dates, where "50009 12/04/97" is selected. The "Series:" column shows "1" selected. The "Image:" column lists images 1 through 8. Below the columns is a "Raw Data Selected:" section with "Yes" and "No" radio buttons, where "No" is selected. At the bottom is a "Comments from User:" text area and a row of buttons: "Save", "Send", "File...", and "Exit". A status bar at the very bottom reads "Loading data completed. Select those you want to save".

Exam:	Series:	Image:
50009 12/04/97	1	1
16 12/04/97		2
15 12/04/97		3
14 12/04/97		4
13 12/04/97		5
12 12/04/97		6
11 12/04/97		7
10 12/04/97		8
9 12/04/97		
8 12/04/97		
7 12/04/97		
50005 12/03/97		
50004 12/03/97		
6 12/03/97		
5 12/03/97		
4 12/03/97		
3 12/03/97		
50003 12/02/97		

Raw Data Selected:
☐ Yes ☒ No


Comments from User:

Save Send File... Exit

Loading data completed. Select those you want to save

EXAM AND SERIES SELECTION
ILLUSTRATION 1-2

Verify that the selected file(s) are correct and select **[Send]**. File(s) will be sent to the ASC. Select **[Dismiss]**, then **[Exit]** to return to normal operation. (see Illustration 1-3) When the transfer to the ASC is successful, the file will automatically be deleted. However, if the files are not sent or an error occurs during transfer, no more than 10 image snap files (saved as .tar files) can be saved. A dialog box will be displayed if more than 10 files exist or if the memory allocation is insufficient to save the selected file. Select the **[File ...]** button to display the files that exist. The files can then be re-sent or deleted as necessary.

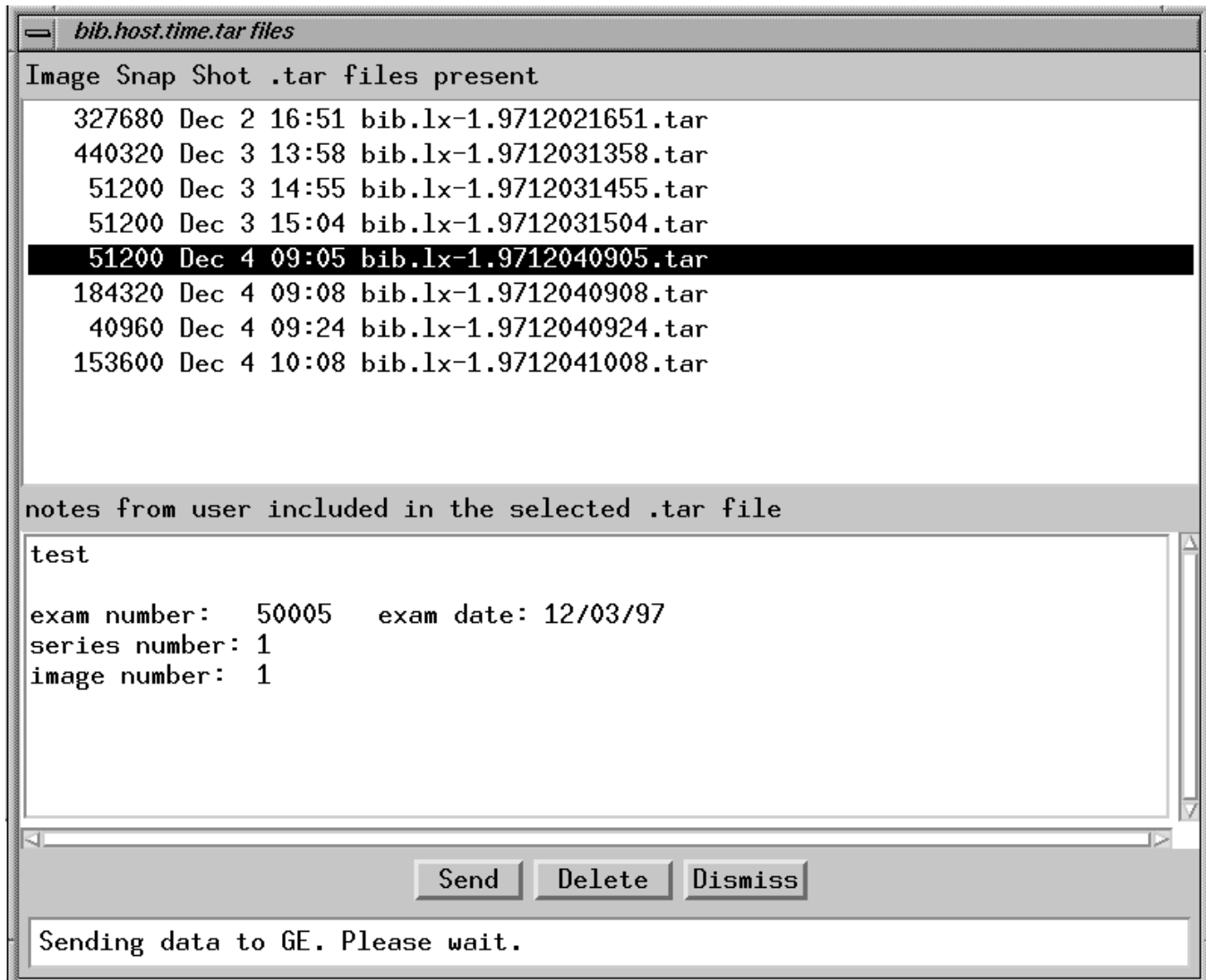

SENDING FILE(S)
ILLUSTRATION 1-3

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
A	12/4/97	W. Pasciak	Initial Release