


TABLE OF CONTENT

TABLE OF CONTENT	1
1- Patient Transport Component Replacement	2
2- Patient Transport Side Cover Removal	3
3- Cradle Assembly Removal.....	4
4- Patient Table Hydraulic Cylinder	5
4-1 Remove Hydraulic Cylinder.....	5
4-2 Install Replacement Hydraulic Cylinder	8
4-3 Final Calibration and Checks	11
5- Patient Table Hydraulic Pump Assembly	14
5-1 Remove Hydraulic Pump Assembly.....	14
5-2 Install Replacement Hydraulic Pump Assembly	20
6- Patient Table Casters	24
7- Dock Hook.....	26
8- Table Up Limit Cable.....	27
9- Secondary Cradle Latch Actuating Cable	30
10- Cradle Interlock Cable.....	34
11- Main Cradle Latch Cable.....	39
12- Hydraulic Filter Replacement	47
12-1 Remove Hydraulic Filter.....	47
12-2 Install Replacement Hydraulic Filter.....	49
13- Bleeding Patient Transport Hydraulics	50
14- Signa Table Arm Board Bumper Field Repair/Replacement Procedure	53
1- Tools and Materials.....	53
2- Repair/Replacement Procedure.....	53
REVISION HISTORY.....	55

Description - Procedures for replacing hardware on the Lightweight Patient Transport.

1- PATIENT TRANSPORT COMPONENT REPLACEMENT

Most replacement procedures covering components of the patient transport require the removal of the transport side covers. Refer to Section 2- Patient Transport Side Cover Removal. Some replacement procedures also require the removal of the cradle. Refer to Section 3- Cradle Assembly Removal. See Illustration L2653A for location of major components of the patient transport.

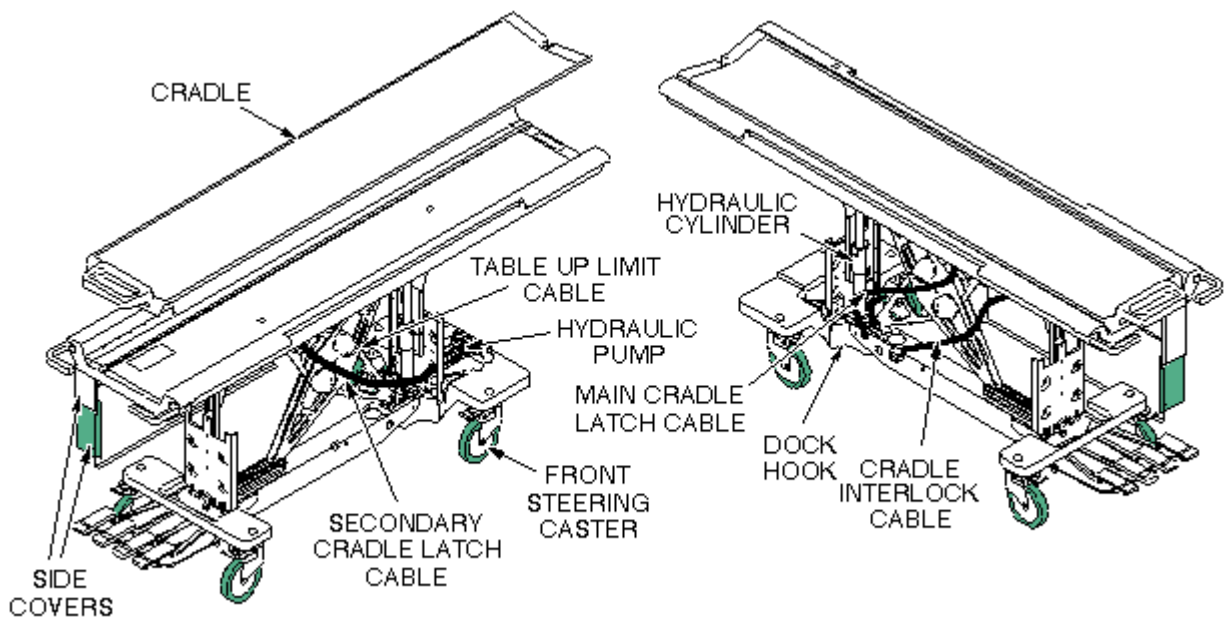


L2653A

LIGHT WEIGHT PATIENT TRANSPORT
ILLUSTRATION L2653A

2- PATIENT TRANSPORT SIDE COVER REMOVAL

See Illustration L2654A for location of upper and lower side covers on patient transport and attaching hardware.


1. Undock patient transport, and move it out of exam room.
2. Remove both outer covers by removing mounting screws (see Illustration L2654A)
3. Remove both inner covers by removing mounting screws.

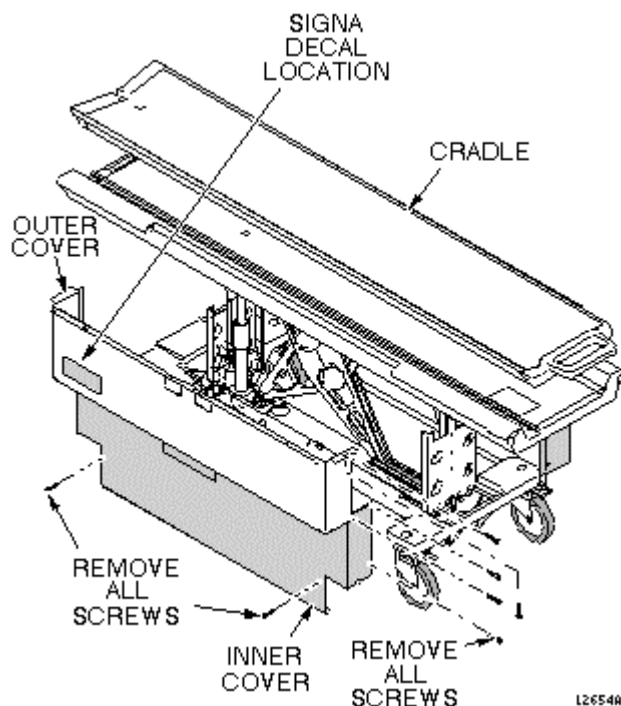

REMOVAL OF PATIENT TABLE COVERS
ILLUSTRATION L2654A

3- CRADLE ASSEMBLY REMOVAL

1. Undock patient transport.
2. Press cradle latch block at front end of transport with a blunt object (such as a pencil eraser), and simultaneously push in actuator rod at front of cradle (see Illustration L2655A)
3. While holding these actuators, slide cradle forward on patient transport approximately three inches. The actuators can now be released.
4. Push in right actuator in base at front end of transport, and slide cradle the rest of the way off of the patient transport (see Illustration L2655A).

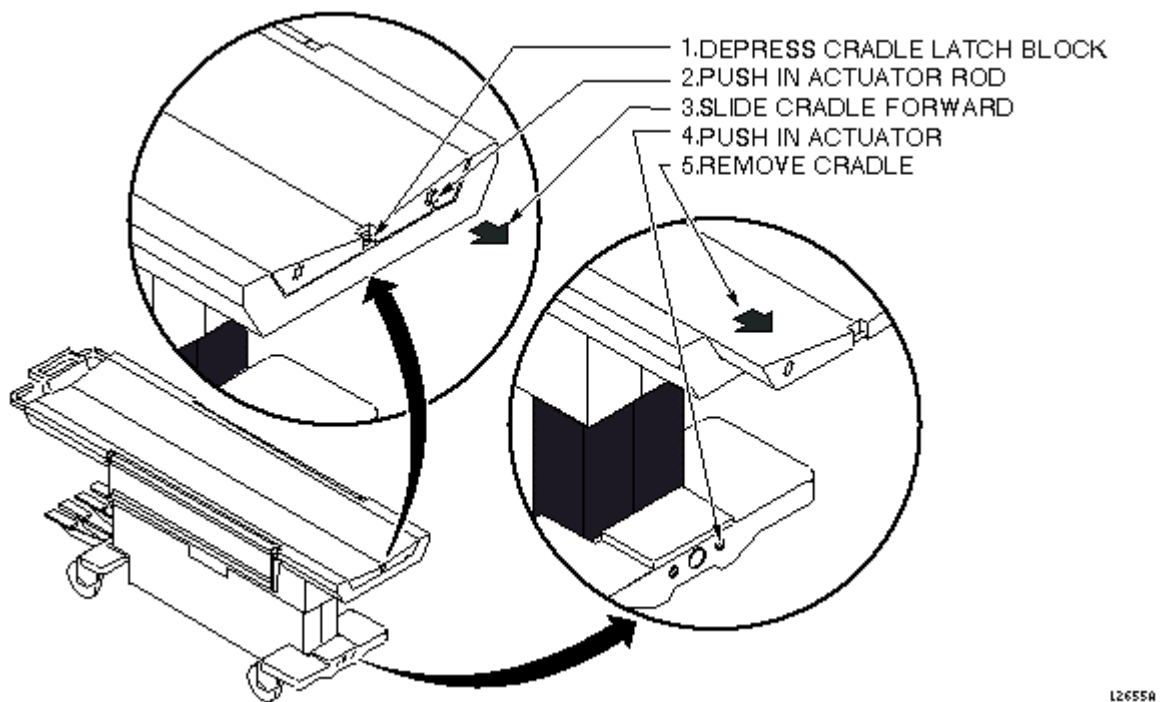

L2655A

PATIENT TRANSPORT CCRADLE REMOVAL
ILLUTRATION L2655A

4- PATIENT TABLE HYDRAULIC CYLINDER


4-1 Remove Hydraulic Cylinder

1. Undock patient transport, and move it out of exam room.
2. Remove both outer covers by removing mounting screws (see Illustration L2654A).
3. Remove both inner covers by removing mounting screws.
4. Remove cradle assembly (see Illustration L2655A).
5. Raise table to *up* position.

WARNING!

PERSONAL INJURY POSSIBILITY! WHEN PERFORMING MAINTENANCE ON THE TABLE IN ITS *UP* POSITION, ENSURE THAT THE LOCK BAR IS SAFELY IN PLACE. SEVERE PERSONAL INJURY CAN OCCUR IF TABLE IS ACCIDENTALLY LOWERED DURING MAINTENANCE PERFORMANCE.


6. Unhook lock bar spring, and lower bar to horizontal position, as shown in Illustration L2656A, then lower table until lock bar fully supports table top.

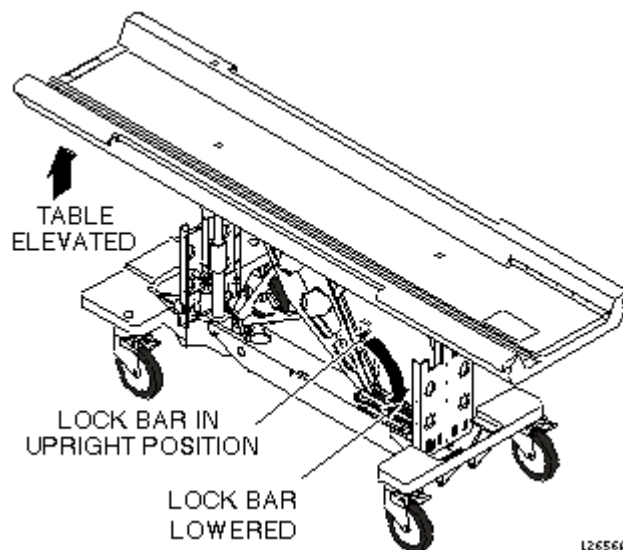

**PATIENT TABLE LOCK BAR
ILLUSTRATION 2656A**

7. Remove plug from patient transport (see Illustration L2657A).

8. Remove retaining ring from end of adjustment screw.
9. Remove retaining ring, washer, and Up Limit cable attachment pin. Then remove tubing and washer. All are in Illustration L2657A.

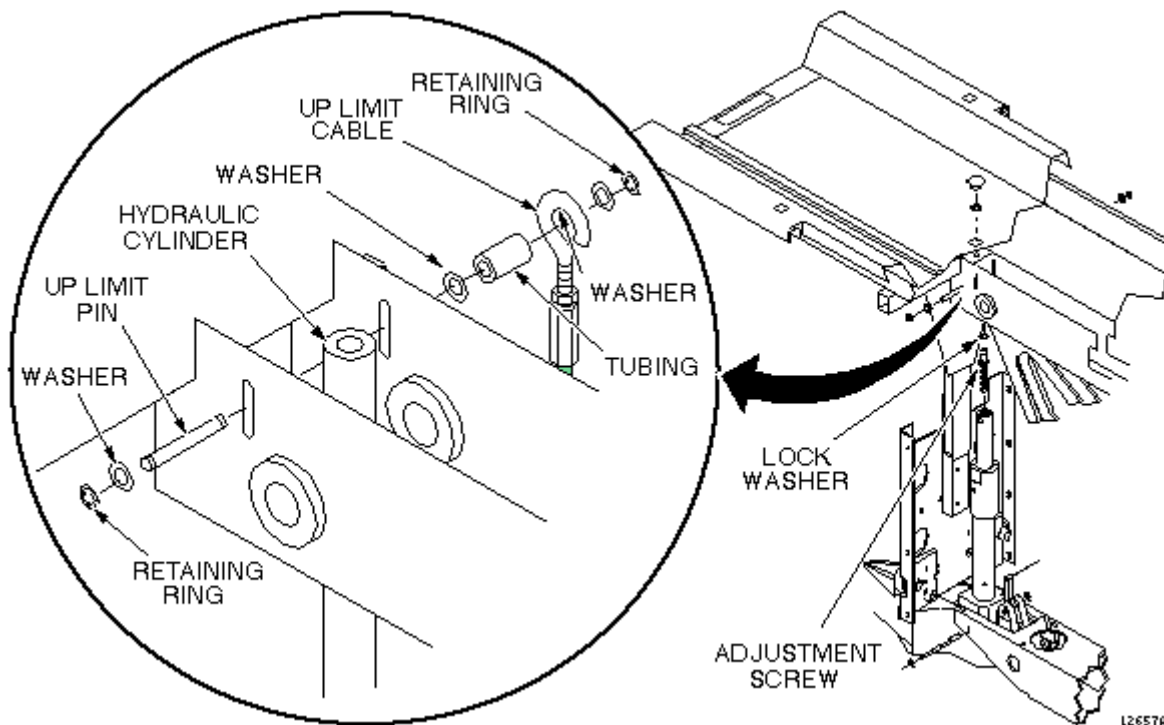

PATIENT TABLE HYDRAULIC CYLINDER REPLACEMENT
ILLUSTRATION L2657A

10. Remove Up Limit cable attachment pin that extends through hydraulic cylinder.


Note

The hydraulic cylinder cannot be pressed unless the Down valve is opened. See Illustration L2659A.

=

11. Manually press hydraulic cylinder down until it is fully retracted.

12. Disconnect two hydraulic lines from base of hydraulic cylinder, and fasten upward and aside to avoid excess spillage. Immediately fold and clamp the translucent (yellow) hose to prevent oil from draining from the reservoir (see Illustration L2660A).
13. Remove retaining ring from either end of cylinder mounting pin at base of cylinder assembly (see Illustration L2660A).

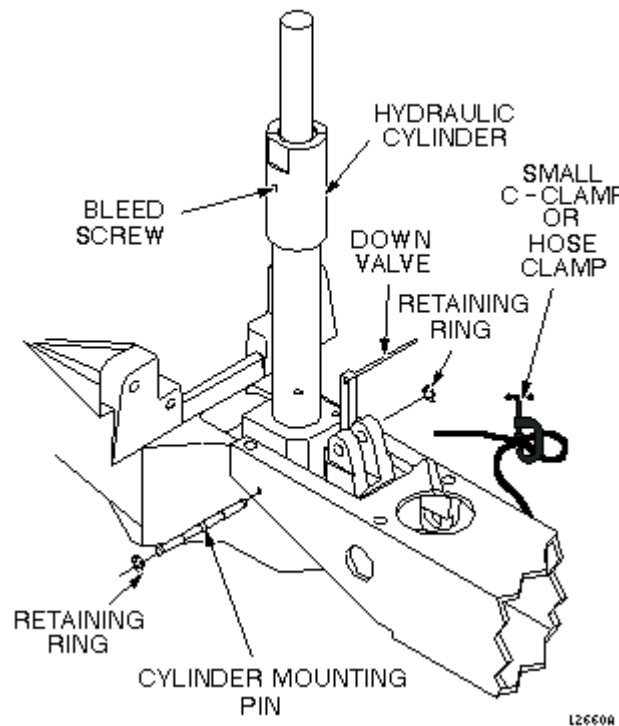

REMOVAL OF HYDRAULIC CYLINDER
ILLUSTRATION L2660A

14. Remove cylinder mounting pin from base of hydraulic cylinder.


15. Carefully move linkage to free actuation pin on Down valve (see Illustration L2661A).
16. Remove hydraulic cylinder from base of table assembly.

4-2 Install Replacement Hydraulic Cylinder

1. Prepare replacement cylinder for installation by pressing cylinder to its fully retracted position.

Note


It may be necessary to move linkage to fit actuation pin of Down valve into groove of linkage, as shown in Illustration L2661A.

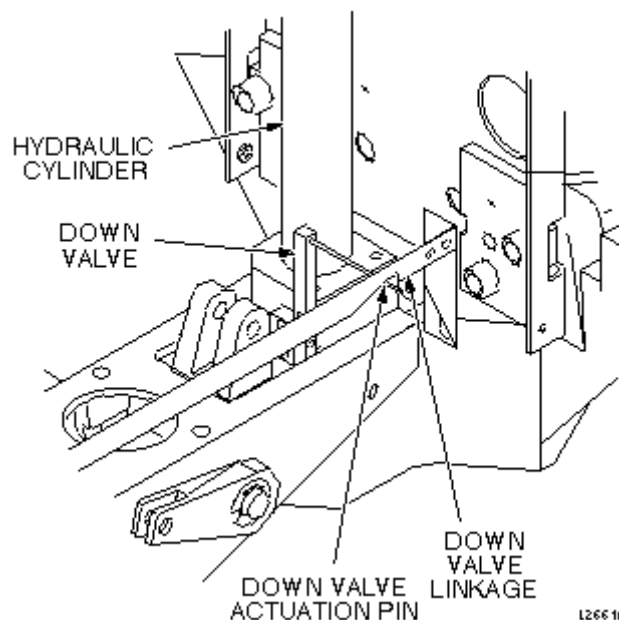

DOWN VALVE ACTUATION PIN
ILLUSTRATION L2661A


2. Carefully set cylinder in base of table, taking care not to bend actuation pin on Down valve.
3. Insert cylinder mounting pin securing base of hydraulic cylinder to frame of table (see Illustration L2660A).
4. Insert retaining ring.
5. Connect two hydraulic lines to base of hydraulic cylinder.

6. Insert adjustment screw into top of hydraulic cylinder (see Illustration L2657A).
7. Insert lock washer on adjustment screw.
8. Position hydraulic cylinder so that top of adjustment screw is aligned with hole in table (see Illustration L2657A).
9. Pump up cylinder until table goes up slightly.
10. Insert retaining ring at top of adjustment screw to secure cylinder to table (see Illustration L2657A).
11. Insert Up Limit cable attachment pin through hole in hydraulic cylinder (see Illustration L2657A). Do not attach the cable at this time.

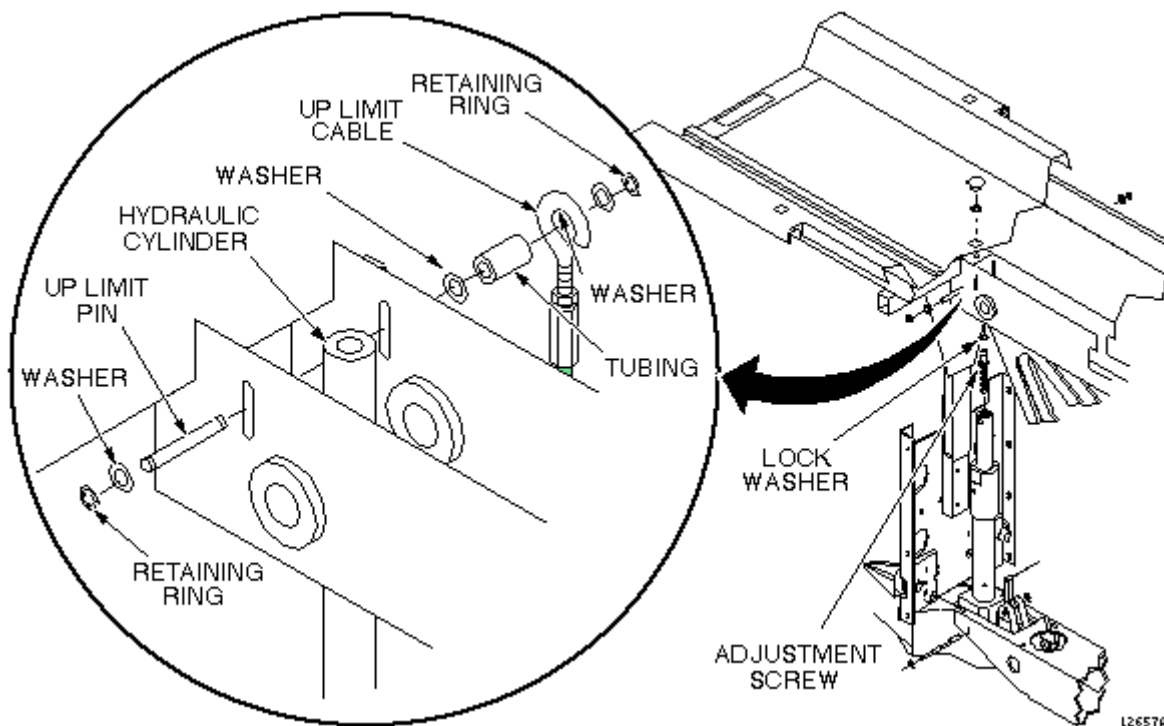


PATIENT TABLE HYDRAULIC CYLINDER REPLACEMENT
ILLUSTRATION L2657A

12. Pump table up to the top of its travel.
13. Attach washer and tubing to Up Limit cable attachment pin.
14. Attach cable to end of Up Limit pin, and secure with washer and retaining ring. Adjust length if necessary to install, and to remove any slack.
15. Ensure that lock bar is at its vertical (Unlocked) position as shown in Illustration L2656A, with the lock bar spring secured.

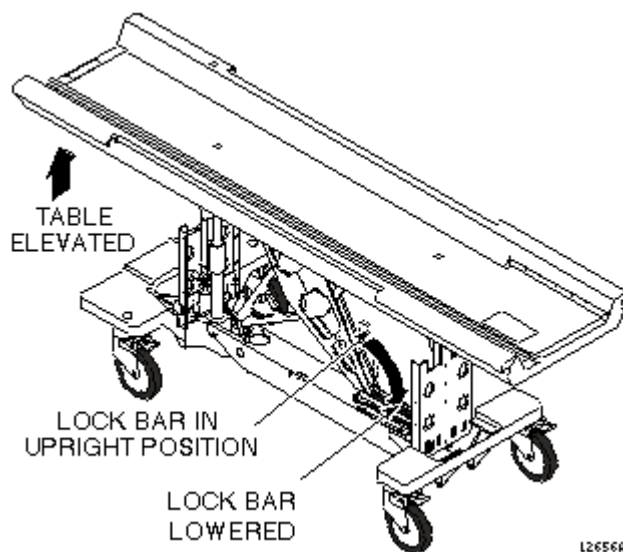

PATIENT TABLE LOCK BAR
ILLUSTRATION 2656A

16. Partially open bleed screw on hydraulic cylinder until no air is visible (see Illustration L2659A).

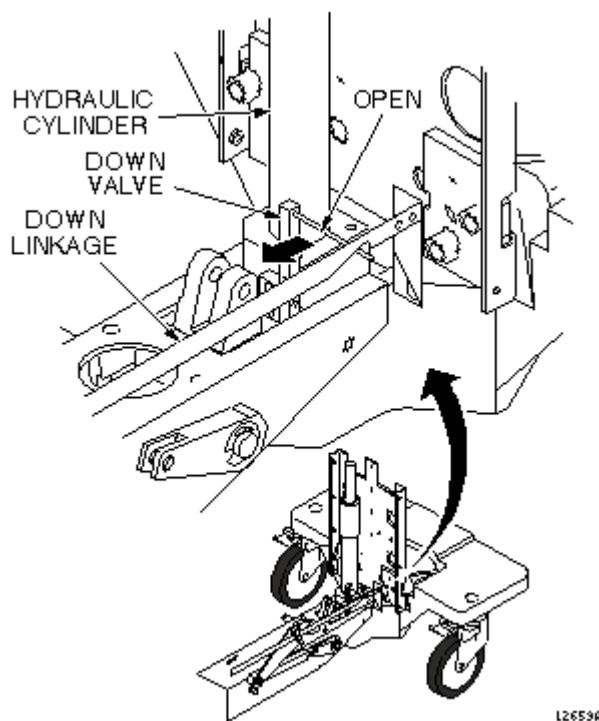


HYDRAULIC CYLINDER DOWN VALVE
ILLUSTRATION L2659A

17. Open Down valve and lower table to its lowest position.

18. Hold Down valve open and exercise Up Pump foot pedal. This bleeds the hydraulic system by shunting air to the hydraulic reservoir.
19. Bleed the cylinder again by partially opening bleed screw, when table is fully up, until no air is visible.


4-3 Final Calibration and Checks

1. Refer to procedure for Lightweight Patient Transport, Section 5, Table Top Height Adjustment; perform those procedures, then perform the following steps.
2. Fill table reservoir so that fluid is 1/4 to 1/2 inch above *full* line when table is at its lowest position (see Illustration L2593A).

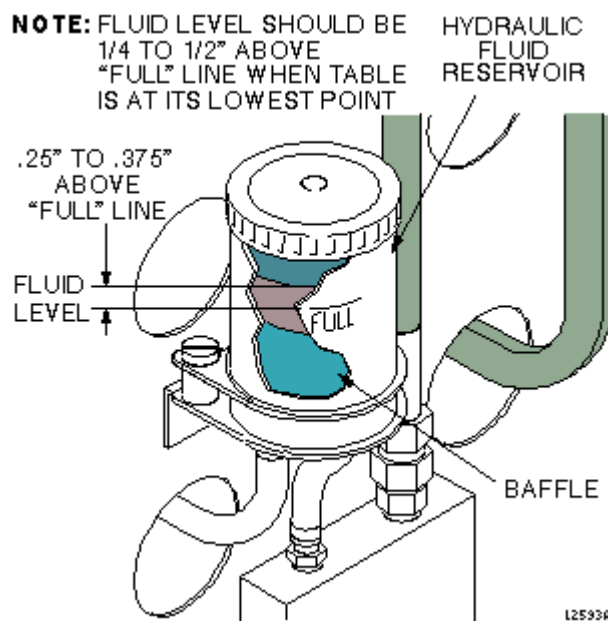


PATIENT TABLE RESERVOIR FILL LEVEL
ILLUSTRATION L2693A

3. Check Down valve linkage adjustment as shown in Illustration L2663A. Obtain desired tolerance between Down valve linkage and bracket by adjusting two set screws on bracket.

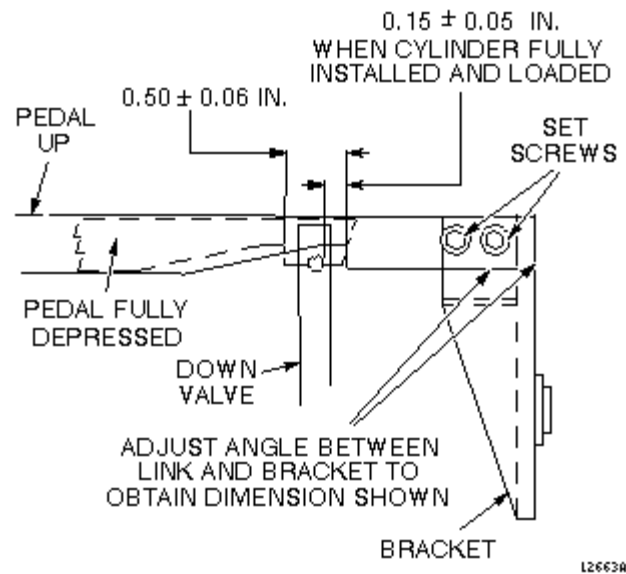


DOWN VALVE LINKAGE ADJUSTMENT
ILLUSTRATION L2663A

4. Replace plug in top of table (see Illustration L2657A).

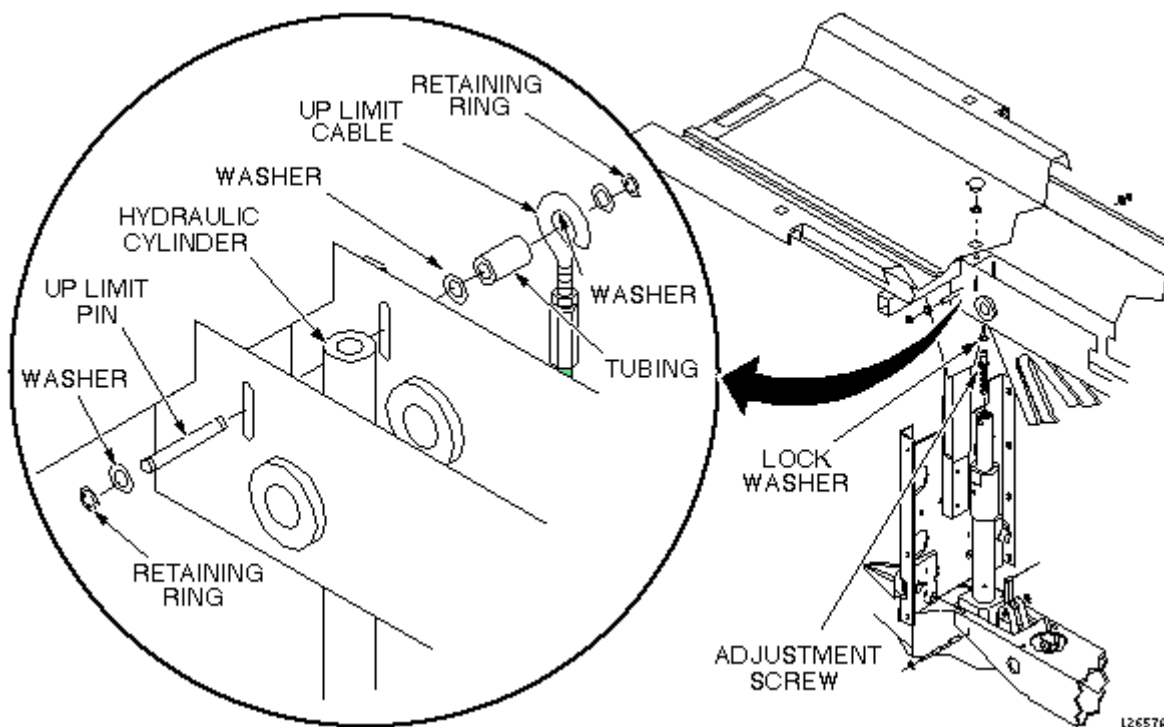


PATIENT TABLE HYDRAULIC CYLINDER REPLACEMENT
ILLUSTRATION L2657A

5. Replace inner and outer covers. Note position of Signa logo on outer covers and position toward front of table (see Illustration L2654A).

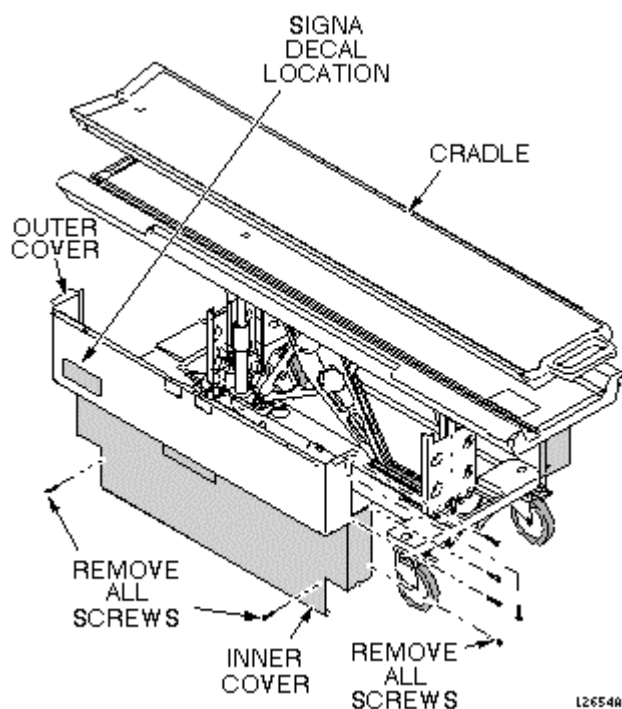

REMOVAL OF PATIENT TABLE COVERS
ILLUSTRATION L2654A

6. Replace cradle assembly (see Illustration L2655A).

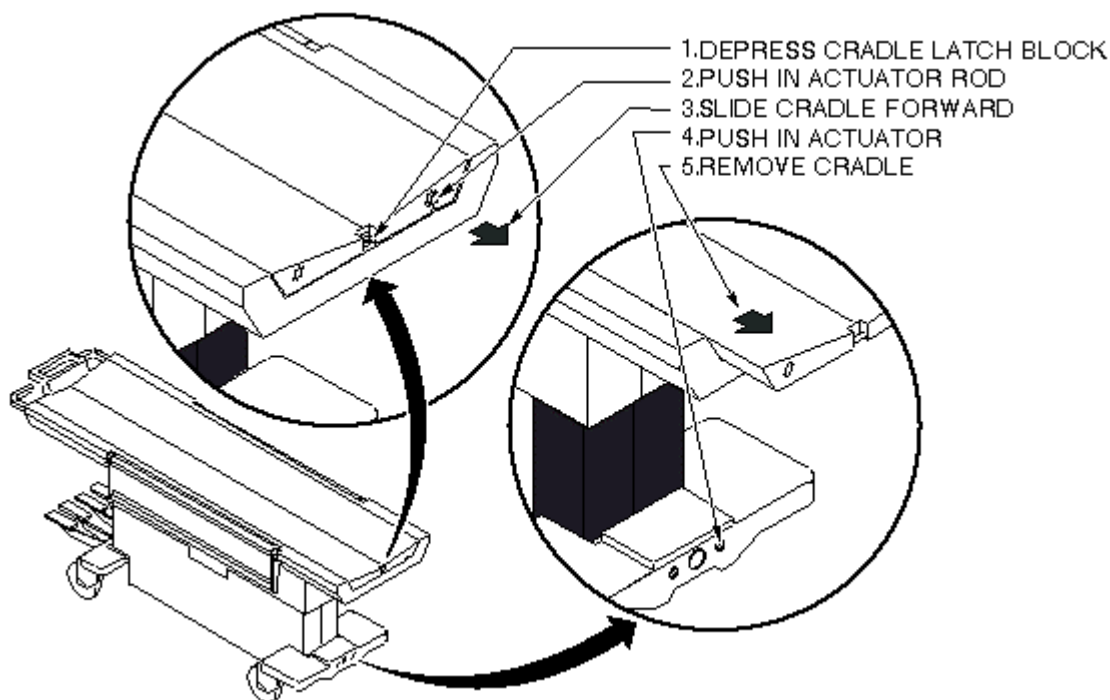

PATIENT TRANSPORT CRADLE REMOVAL
ILLUSTRATION L2655A


5- PATIENT TABLE HYDRAULIC PUMP ASSEMBLY

5-1 Remove Hydraulic Pump Assembly

Note

Have sufficient wipes and a container available to catch hydraulic fluid.

1. Undock patient transport, and move it out of exam room.
2. Remove both outer covers by removing mounting screws (see Illustration L2654A).

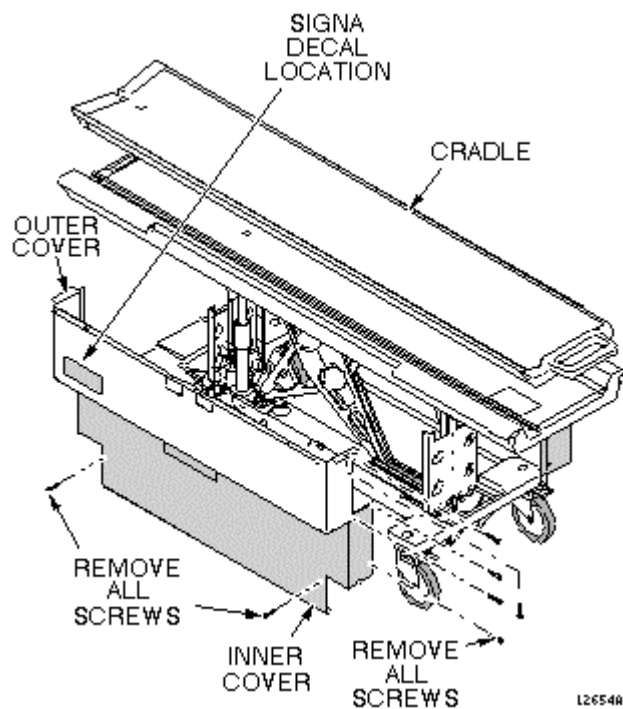


REMOVAL OF PATIENT TABLE COVERS
ILLUSTRATION L2654A

3. Remove inner covers by removing mounting screws.
4. Raise table to *up* position.

WARNING!


PERSONAL INJURY POSSIBILITY! WHEN PERFORMING MAINTENANCE ON THE TABLE IN ITS *UP* POSITION, ENSURE THAT THE LOCK BAR IS SAFELY IN PLACE. SEVERE PERSONAL INJURY CAN OCCUR IF TABLE IS ACCIDENTALLY LOWERED DURING MAINTENANCE PERFORMANCE.

5. Unhook lock bar spring and move lock bar to horizontal position, as shown in Illustration L2656A, then lower table until lock bar is supporting table in its full up position.

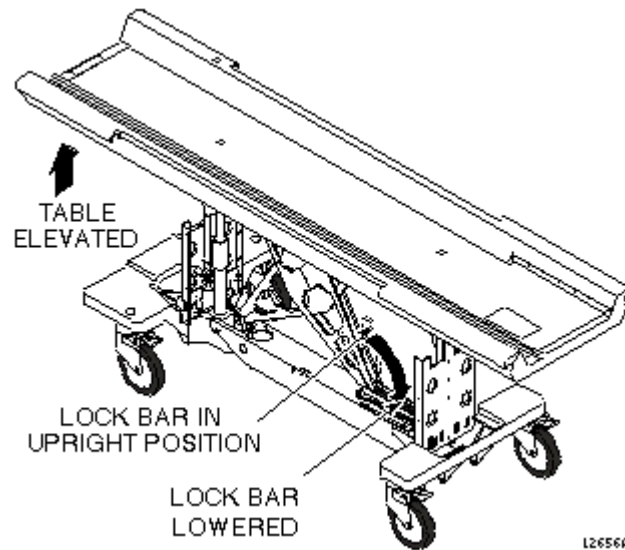


PATIENT TABLE LOCK BAR
ILLUSTRATION L2656A

6. Disconnect the two hydraulic lines to hydraulic pump assembly, and fasten upward and aside to avoid excess spillage. Immediately fold and clamp the translucent (yellow) hose from reservoir to prevent oil from draining from the reservoir (see Illustration L2664A).

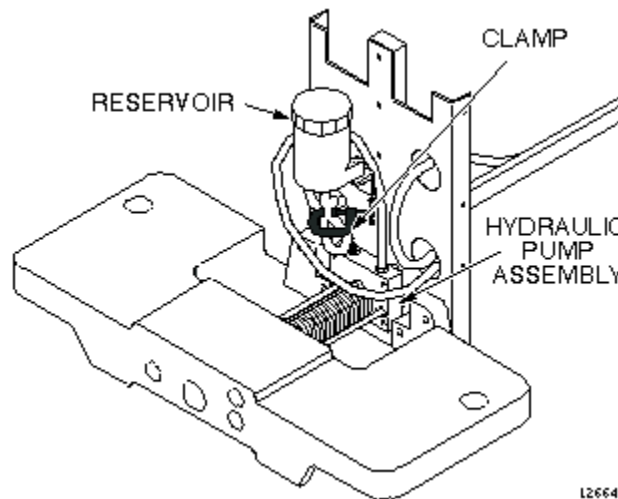

HYDRAULIC PUMP ASSEMBLY
ILLUSTRATION L2664A


7. Loosen two mounting screws attaching upper drip pan to lower drip pan (see Illustration L2665A).
8. Remove upper drip pan (see Illustration L2665A).

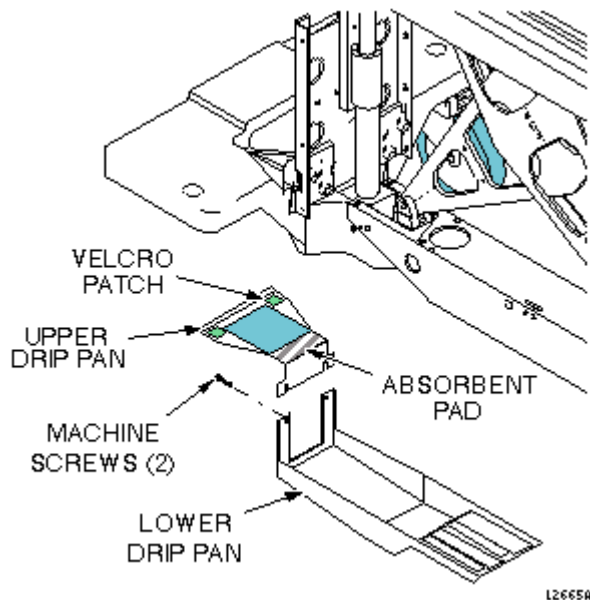

UPPER AND LOWER DRIP PANS
ILLUSTRATION L2665A


CAUTION

Equipment damage possibility. Do not remove spring pins from either end of pump pull rods; this will cause a side loading on the pump, requiring its replacement. The pump pull rods are disconnected from the pump linkage by turning the pump pull rods in a counterclockwise direction after tension has been removed from the pump springs. This is explained in the following procedure. The pump pull rod spring pins are exposed only when the Table Up pedal is pressed (see Illustration L2666A).

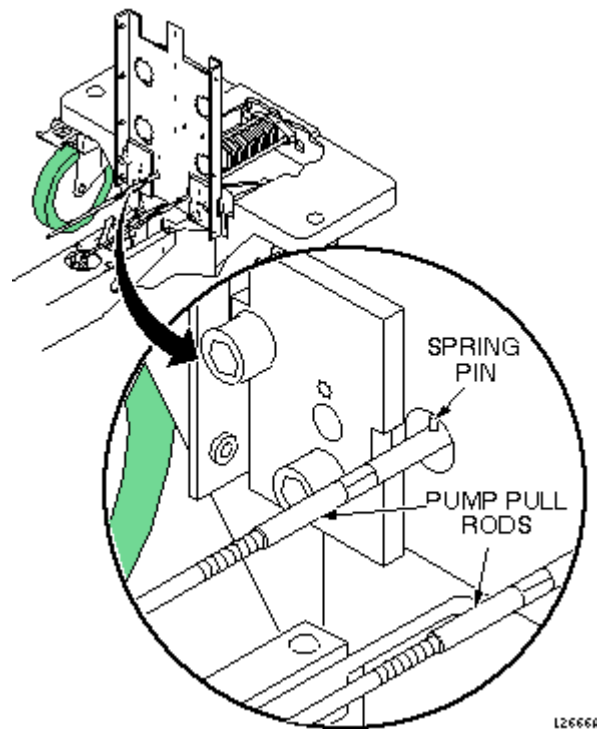


HYDRAULIC PUMP PULL RODS
ILLUSTRATION L666A

9. Compress main pump spring, as shown in Illustration L2667A, to remove tension from pull rods. Block the pump Up pedal fully up to take its weight off of the pull rods.
10. Push the yoke back and unscrew the pull rods from the pump linkage (see Illustration L2666A).

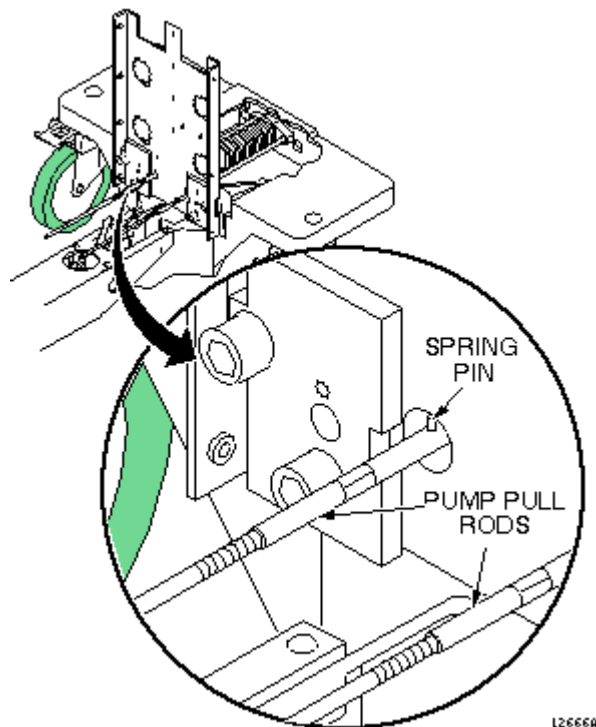


HYDRAULIC PUMP PULL RODS
ILLUSTRATION L666A

11. Remove c-clamps used to compress main spring, as shown in Illustration L2667A.

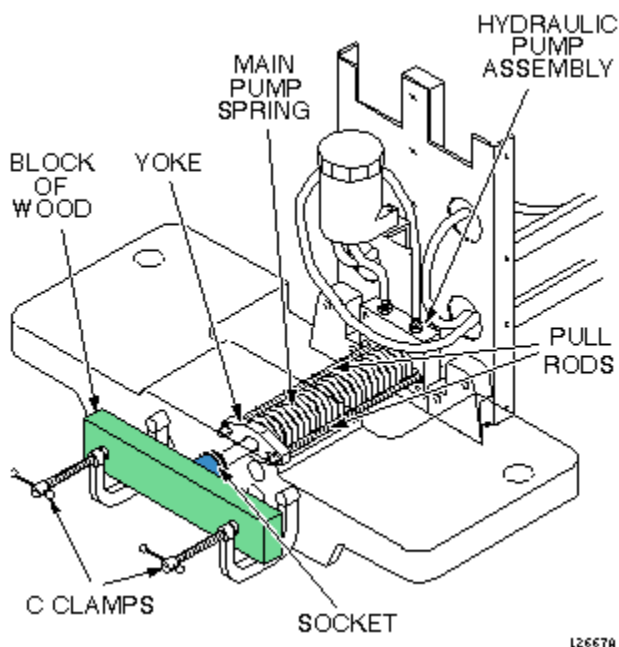

COMPRESS HYDRAULIC PUMP MAIN SPRING
ILLUSTRATION L2667A

12. Remove machine screw securing bumper to the pump piston (see Illustration L2668A).

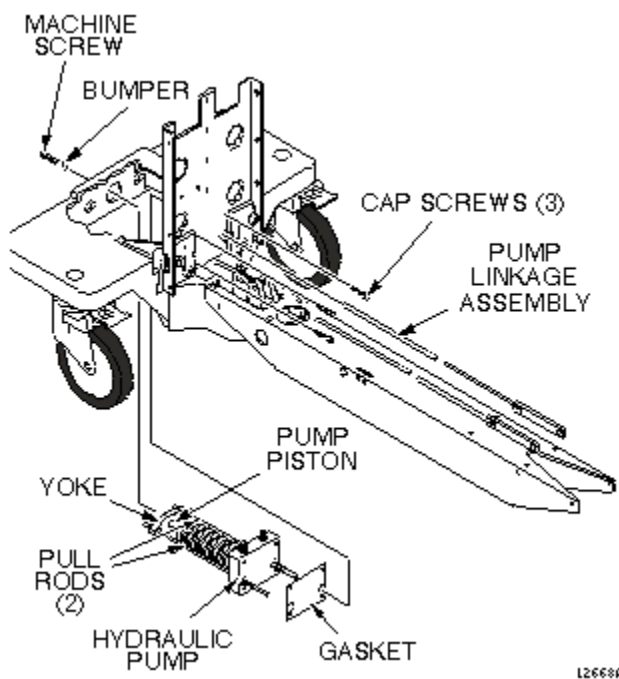


HYDRAULIC PUMP REPLACEMENT
ILLUSTRATION L2668A

13. Remove three mounting cap screws that secure pump assembly to front table support.


14. Remove hydraulic pump assembly and gasket.
15. If absorbent pad attached to upper drip pan appears to be saturated with fluid, remove it and squeeze out fluid.

5-2 Install Replacement Hydraulic Pump Assembly

1. Prepare replacement pump for installation by applying Locktite 242 to three mounting screws and to machine screw that secures bumper to pump piston.
2. Set replacement pump assembly on front support of table, making sure that gasket is installed properly, and secure with three mounting screws (see Illustration L2668A).
3. Secure the bumper to the pump piston with machine screw.
4. Reinstall c-clamps on front support, as shown in Illustration L2667A, to compress main pump spring. Screw pump pull rods into pump linkage, alternately tightening each side to avoid side loading the pump. Stop when the rear spring pin on each pull rod is 1/8-inch (3.2 mm) from pump body (flush with back surface of front support) (see Illustration L2669A).

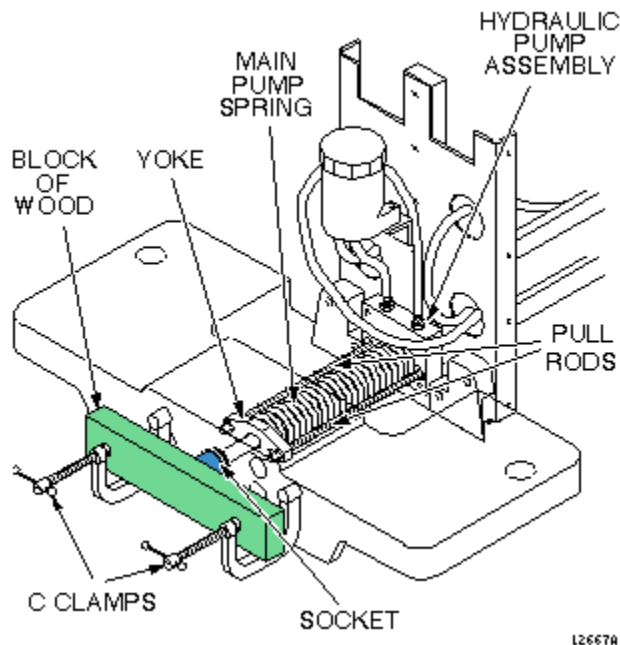


COMPRESS HYDRAULIC PUMP MAIN SPRING
ILLUSTRATION L2667A

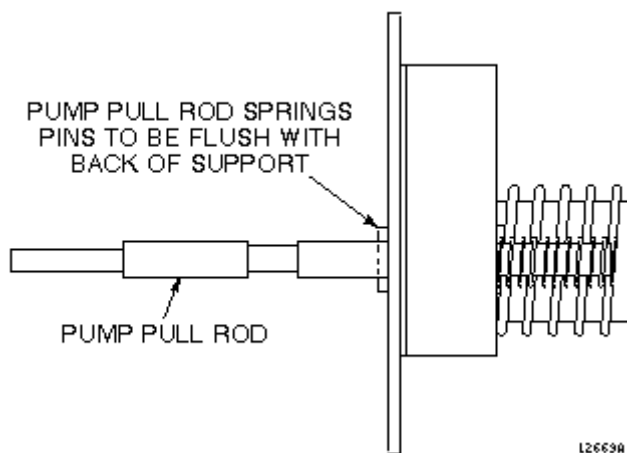


PUMP PULL ROD ADJUSTMENT
ILLUSTRATION L2669A

5. Connect two hydraulic lines to hydraulic pump assembly (see Illustration L2664A).

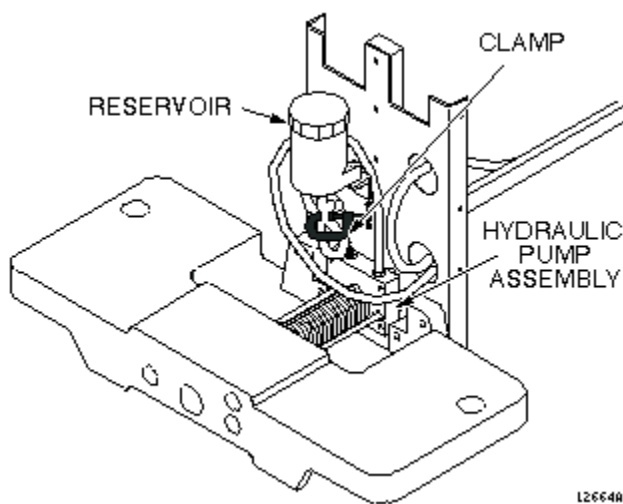

HYDRAULIC PUMP ASSEMBLY
ILLUSTRATION L2664A

6. Attach upper and lower drip pans (see Illustration L2665A).

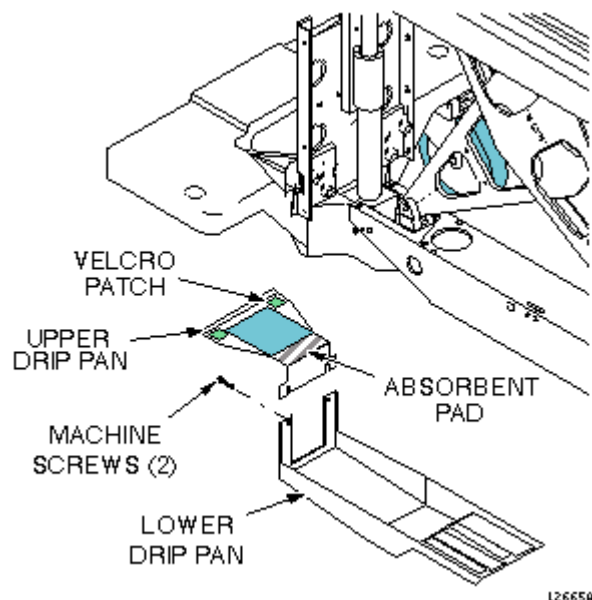


UPPER AND LOWER DRIP PANS
ILLUSTRATION L2665A

7. Block open the Down valve, and pump the Up pedal at least 20 full strokes to purge air from the lines.
8. Close Down valve and raise table.
9. Retract lock bar and secure with spring, then lower table.
10. Fill table reservoir so that fluid is 1/4 to 1/2 inch above *full* line when table is at its lowest position.

Note


If any air gets into the cylinder, open the hydraulic cylinder bleed screw, when the table is raised, until air is removed from the cylinder. Very little fluid loss should occur during a normal bleed procedure.

11. Verify that pump pull rods are adjusted equally by comparing position of rear spring pins on pump pull rods (see Illustration L2669A).

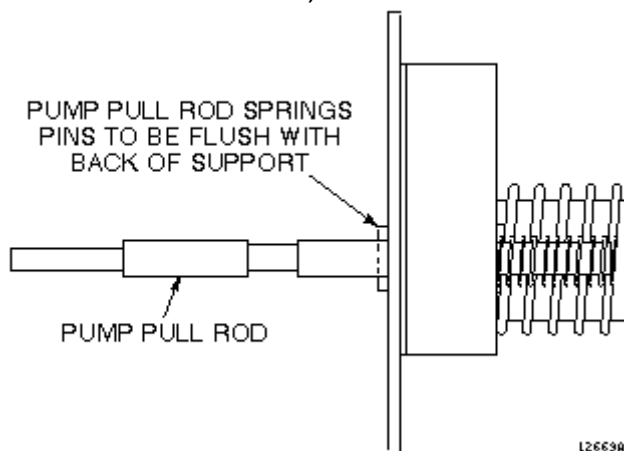

PUMP PULL ROD ADJUSTMENT
ILLUSTRATION L2669A

12. Replace inner and outer covers. Note position of Signa logo on outer covers and position toward front of table (see Illustration L2654A).

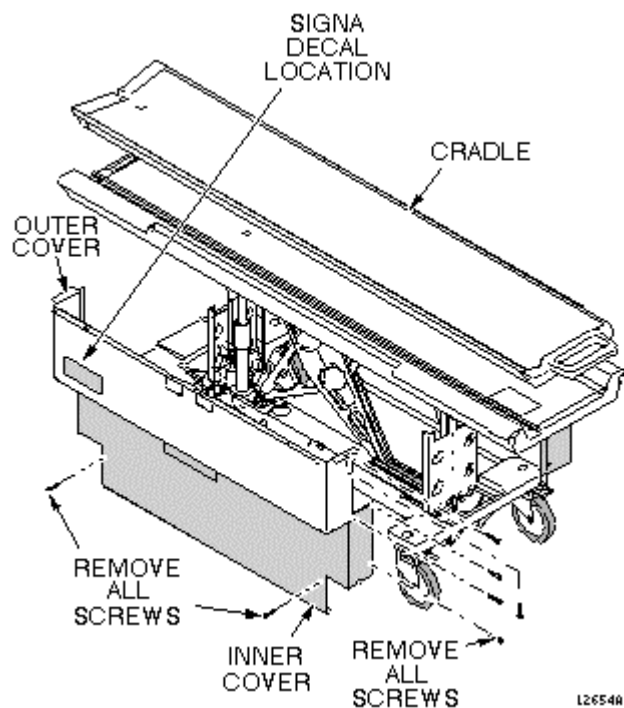

REMOVAL OF PATIENT TABLE COVERS
ILLUSTRATION L2654A

6- PATIENT TABLE CASTERS


There are two types of casters on the patient transport: a straight-line tracking (steering) caster that is mounted on the right front side of the patient transport, or three lock type casters (see Illustration L2670A).

Field-supplied Materials

- Synthetic Lubricating Grease 46-256047P1
- Block of wood to support table when casters are being replaced.

Procedure

1. Unlock patient transport, and move from exam room.
2. Prepare block of wood to support table when casters are replaced.
3. Remove rubber cap from the caster to be replaced (see Illustration L2670A).

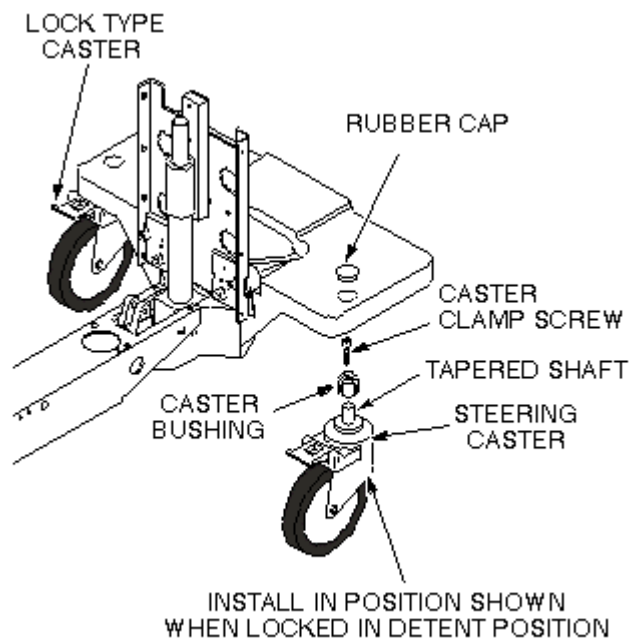

L2670A

CASTER WHEEL ASSEMBLY
ILLUSTRATION L2670A


4. Remove caster clamp screw (center screw) from the caster to be replaced.
5. Repeat steps 3 and 4 on any one of the other casters, because two caster clamp screws are required to force the tapered caster shaft out of the caster bushing on the defective caster.

6. Return to the caster to be replaced. Thread the two caster clamp screws in the vacant caster release screw holes of the caster bushing. Force the tapered caster shaft out of the caster bushing by tightening two caster clamp screws, equally, a little at a time, until caster taper fit breaks loose and caster falls off.
7. Be sure that bushing is not bell-mouthed at bottom. If bushing is damaged, remove by turning out of the bottom of the casting using a 3/4-inch (19 mm) socket.

Note

If new bushing is required, apply synthetic lubricating grease to inside taper and threads of bushing before installation.

8. Using synthetic lubricating grease, lubricate tapered stem of caster and threads of caster clamp screw.
9. Install new caster into bushing. Straight line tracking (steering) caster should be installed so that it is oriented in position shown when locked in detent position (see Illustration L2670A).

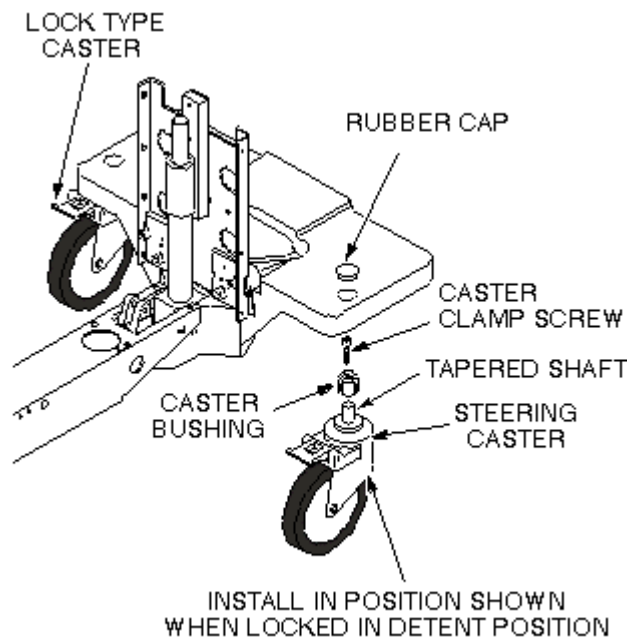

L2670A

CASTER WHEEL ASSEMBLY
ILLUSTRATION L2670A


10. Install caster clamp screw and tighten.
11. Remove support wood blocks.
12. Perform procedure for Lightweight Patient Transport Section 2-3, Checking Caster Height Adjustment.



Equipment damage possibility. Caster adjustments must be done according to section 2-3, Procedure Checking Caster Height Adjustment. Misadjustments leave caster bushing vulnerable to damage, which can result in injury to patient and/or persons near patient transport.

7- DOCK HOOK

1. Undock patient transport, and move from exam room.
2. Loosen set screw under arm assembly (see Illustration L2671A).
3. Remove hex jam nut from the hook and then unscrew hook from the arm assembly (see Illustration L2671A).

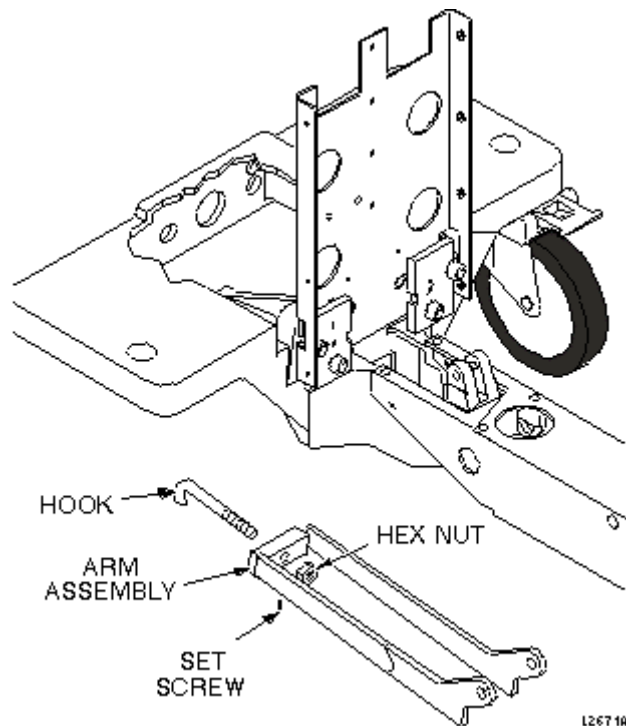

PATIENT TABLE DOCK HOOK
ILLUSTRATION L2671A

4. Install new dock hook, but do not lock down until adjustment is complete.
5. Move patient transport into exam room.

Note

If dock hook appears deformed or stretched after extended use, this indicates that it may be yielding, which can lead to fatigue. If this occurs, replace the hook, rather than tightening it.


WARNING!

EQUIPMENT DAMAGE POSSIBILITY! IF THE DOCK HOOK IS ADJUSTED TOO TIGHTLY, AND AN ATTEMPT IS MADE TO DOCK THE TABLE, SEVERE EQUIPMENT DAMAGE CAN OCCUR. IF IT APPEARS TO TAKE AN EXCESSIVE AMOUNT OF FORCE TO PRESS THE DOCK PEDAL, DO NOT FORCE IT; LOOSEN THE DOCK HOOK FIRST.

6. Perform dock hook adjustment procedures per Set Up And Calibration tab, section 3, Procedure Dock Hook Tension Adjust

8- TABLE UP LIMIT CABLE


1. Undock patient transport, and move from exam room.
2. Remove upper and lower side covers (see Illustration L2654A).

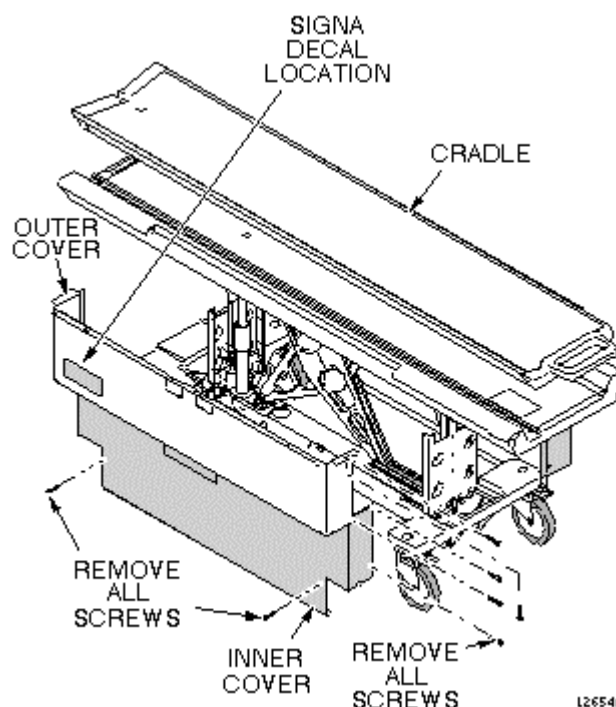

REMOVAL OF PATIENT TABLE COVERS
ILLUSTRATION L2654A

3. Raise table to maximum height.
4. Remove cotter pin and retaining pin securing bottom end of cable to bell crank (see Illustration L2672A).
5. Remove clamp securing cable to bracket (see Illustration L2672A).

6. Remove retaining ring and washer securing end of cable to Up Limit pin (see Illustration L2672A).
7. Install new cable and attach to Up Limit pin per step 6 above.
8. Fasten other end of cable to bell crank, and secure retaining pin with cotter pin (see Illustration L2672A).
9. Attach clamp securing cable to bracket (see Illustration L2672A).

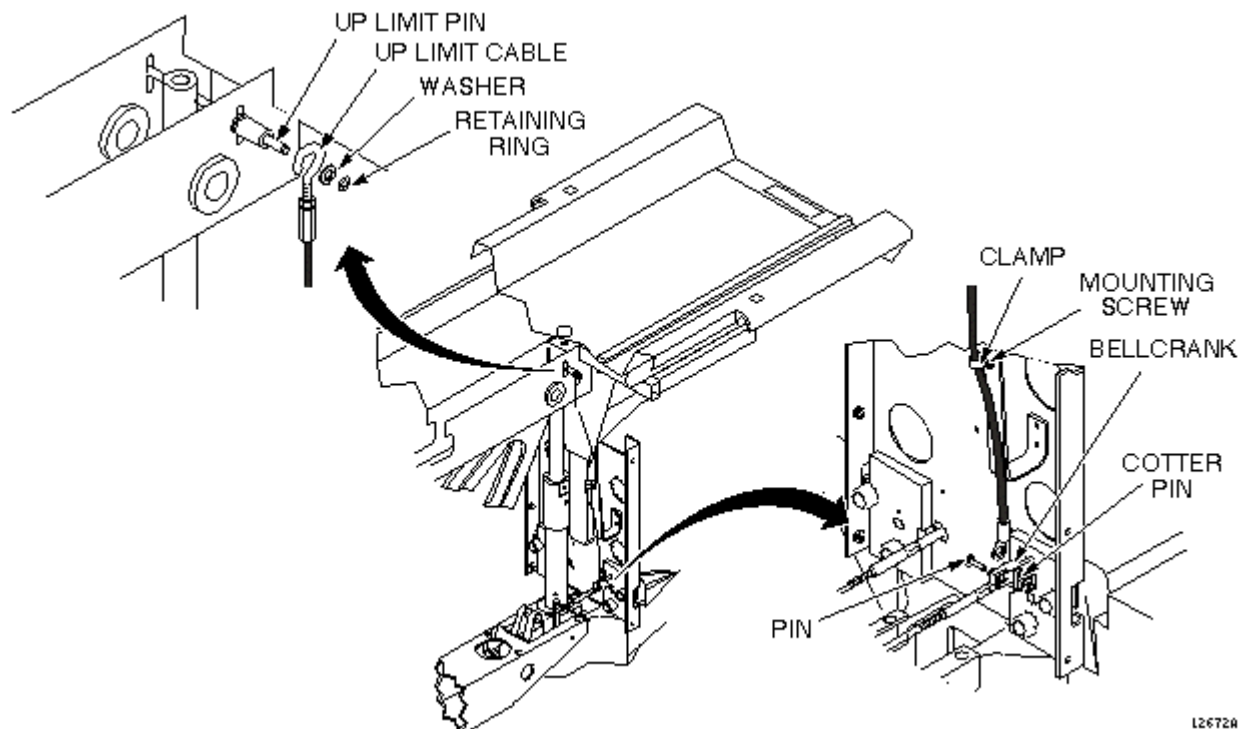


L2672A

TABLE UP LIMIT CABLE REPLACEMENT
ILLUSTRATION L2672A

10. Adjust table Up Limit cable to remove any slack when hydraulic cylinder is fully extended.
11. Check for proper operation of Up Limit switch when table is docked.


12. Replace upper and lower side covers. Note position of Signa logo on outer covers and position toward front of table (see Illustration L2654A).

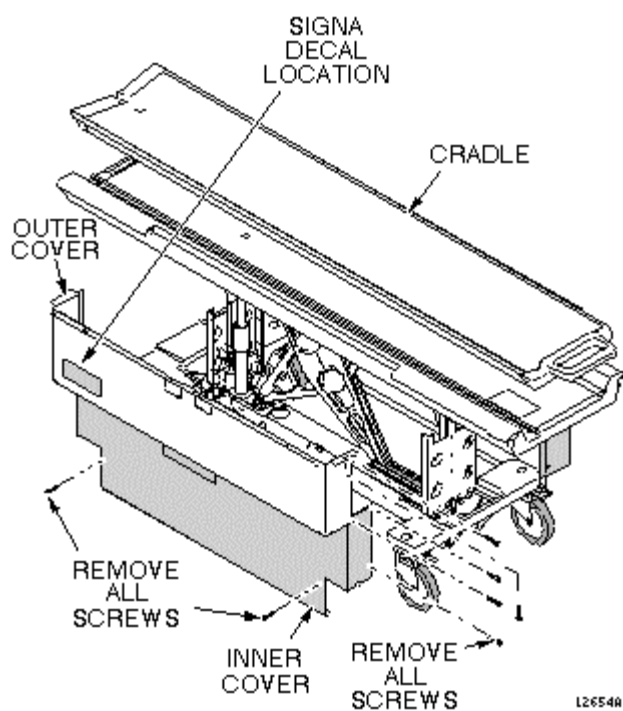

REMOVAL OF PATIENT TABLE COVERS
ILLUSTRATION L2654A


9- SECONDARY CRADLE LATCH ACTUATING CABLE

1. Undock patient transport, and move from exam room (see Illustration L2673A).

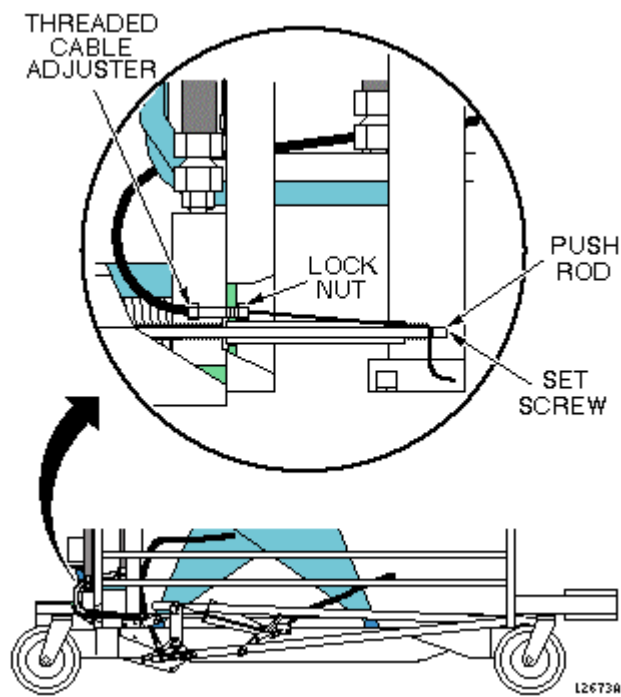


SECONDARY CRADLE LATCH ACTUATING CABLE
ILLUSTRATION L2673A

2. Remove upper and lower side covers (see Illustration L2654A).

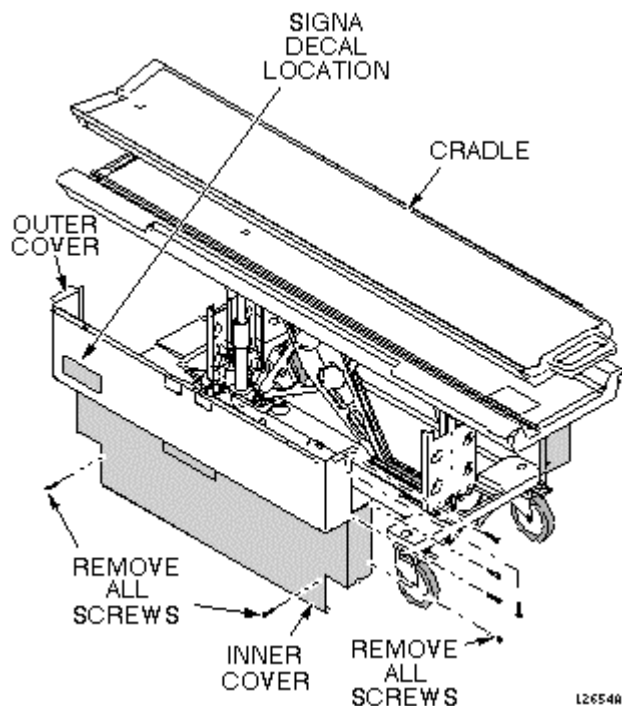


REMOVAL OF PATIENT TABLE COVERS
ILLUSTRATION L2654A

3. Remove cradle (see Illustration L2655A).

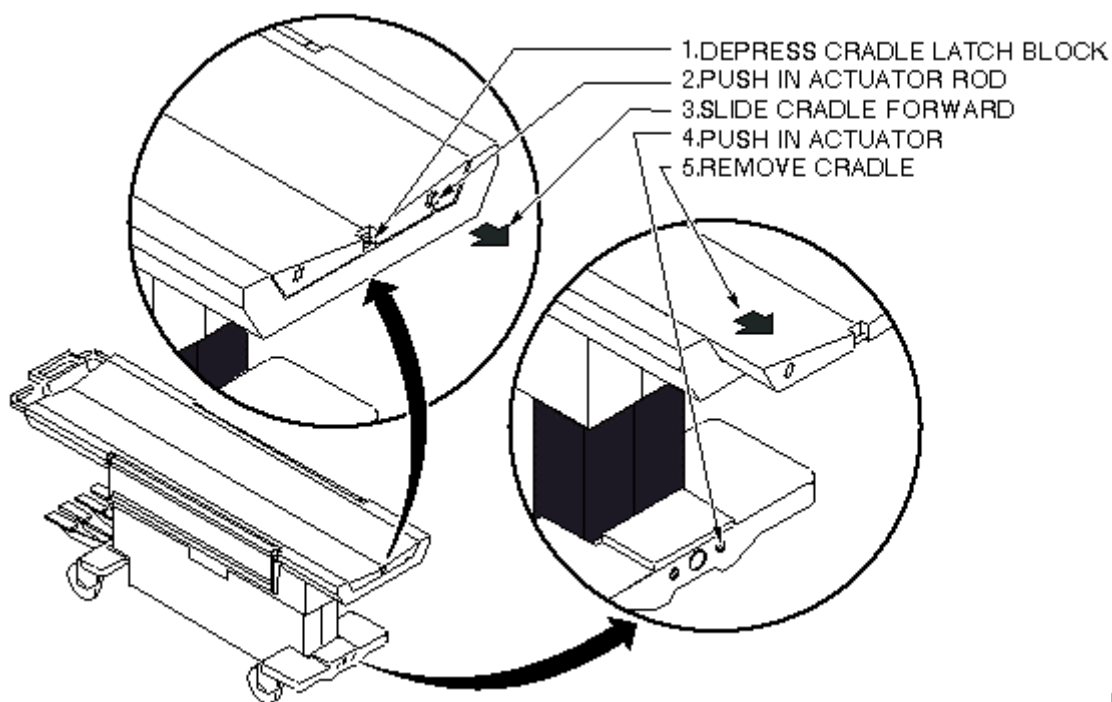


PATIENT TRANSPORT CRADLE REMOVAL
ILLUSTRATION L2655A

4. Unfasten cable end from end of push rod (see Illustration L2673A).

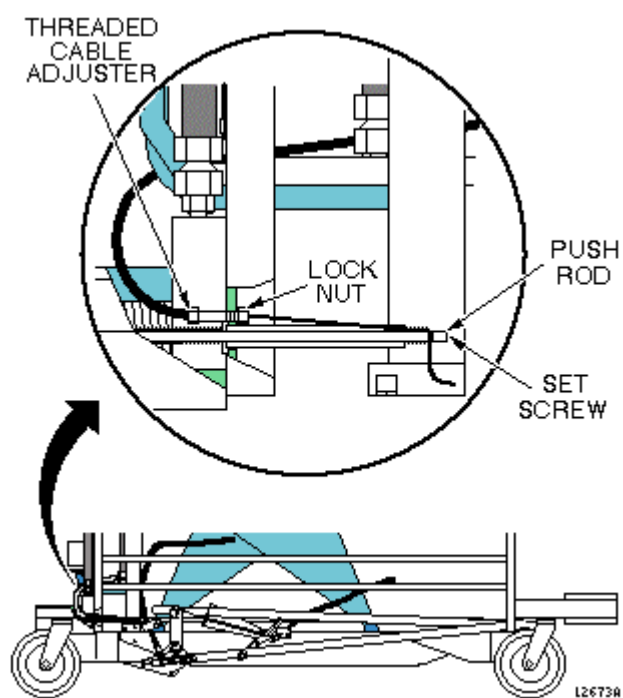

SECONDARY CRADLE LATCH ACTUATING CABLE
ILLUSTRATION L2673A

5. Pull cable out of pin housing (see Illustration L2674A).
6. Insert new cable through cable casing (see Illustration L2674A).
7. Route cable through cable casing (see Illustration L2674A).
8. Loop cable twice through hole in push rod and secure with set screw. Tighten securely to prevent cable from slipping (see Illustration L2674A).
9. Adjust cable so that when push rod is pushed in $\frac{3}{8}$ inch (9.5 mm); pin stop is flush with pin housing. Use threaded cable adjuster for fine adjustment. Tighten lock nut after adjustment is complete, then cut off excess length of cable (see Illustration L2674A).

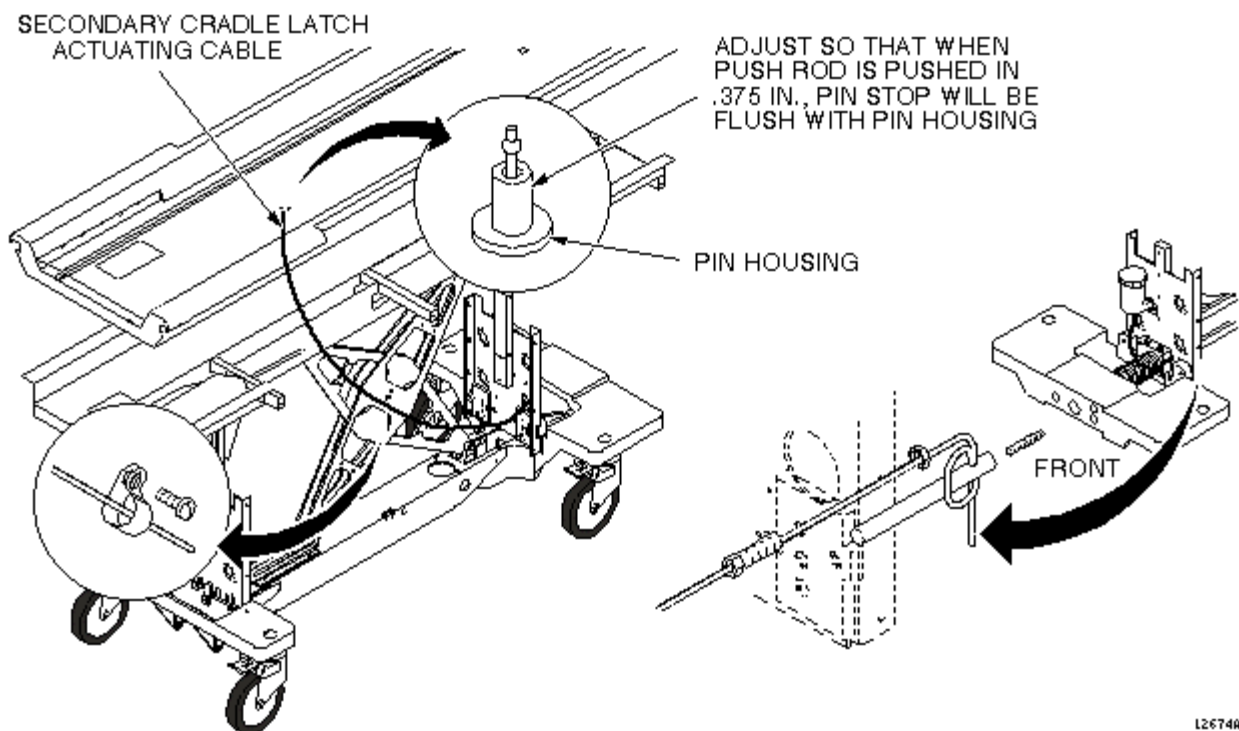


SECONDARY CRADLE LATCH CABLE ROUTING
ILLUSTRATION L2674A

10. Return patient transport to scan room.
11. Dock patient transport.
12. Verify that release pin is flush to $\frac{1}{16}$ inch (1.6 mm) below flush when table is docked.

13. Replace cradle and upper and lower side covers. Note position of Signa logo on outer covers, and position toward front of table (see Illustration L2654A).


10- CRADLE INTERLOCK CABLE

1. Undock patient transport, and move it out of exam room.
2. Remove upper and lower side covers (see Illustration L2654A).

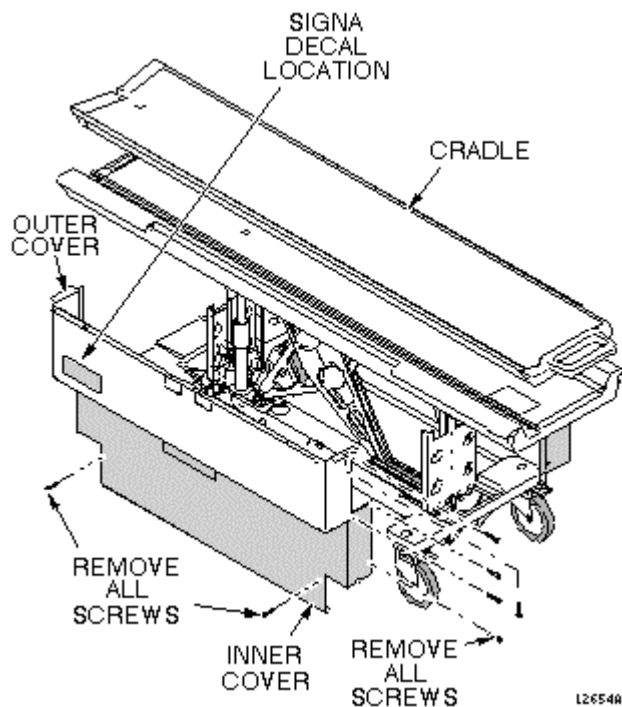

REMOVAL OF PATIENT TABLE COVERS
ILLUSTRATION L2654A


3. Remove small cover below flipper assembly (see Illustration L2675A).

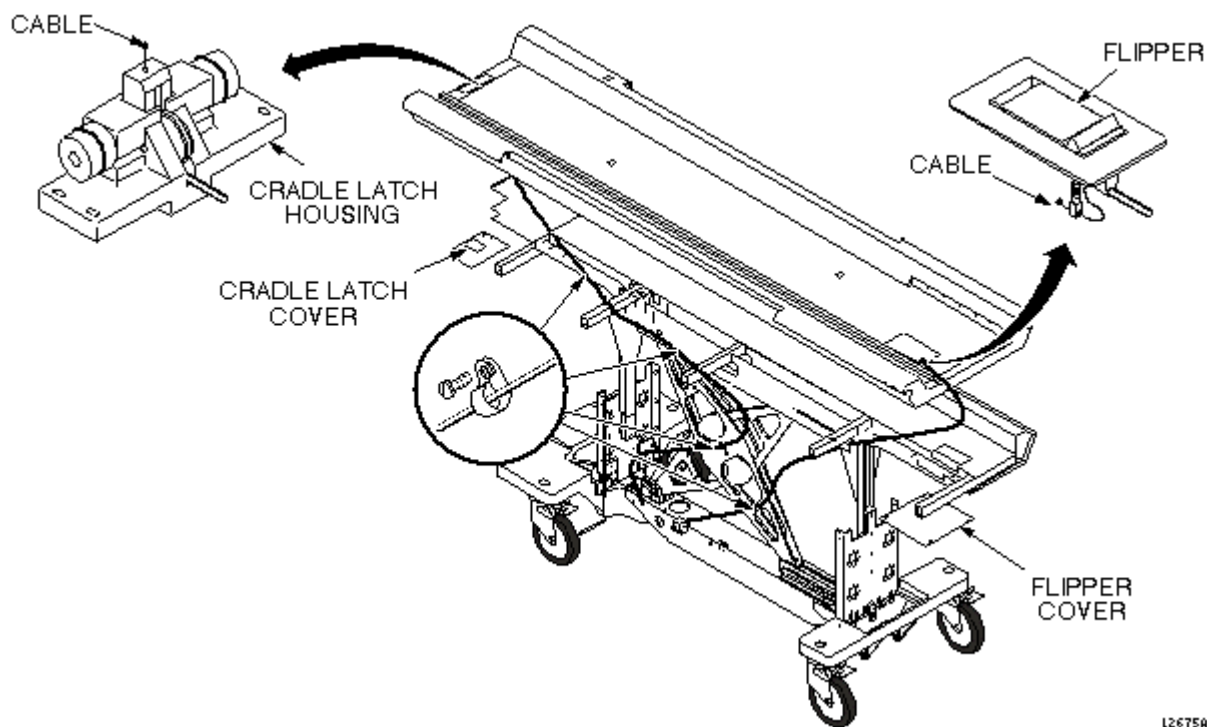

L2675A

INTERLOCK CABLE AND MAIN CRADLE LATCH CABLE ROUTING
ILLUSTRATION L2675A


4. Remove cradle (see Illustration L2655A).

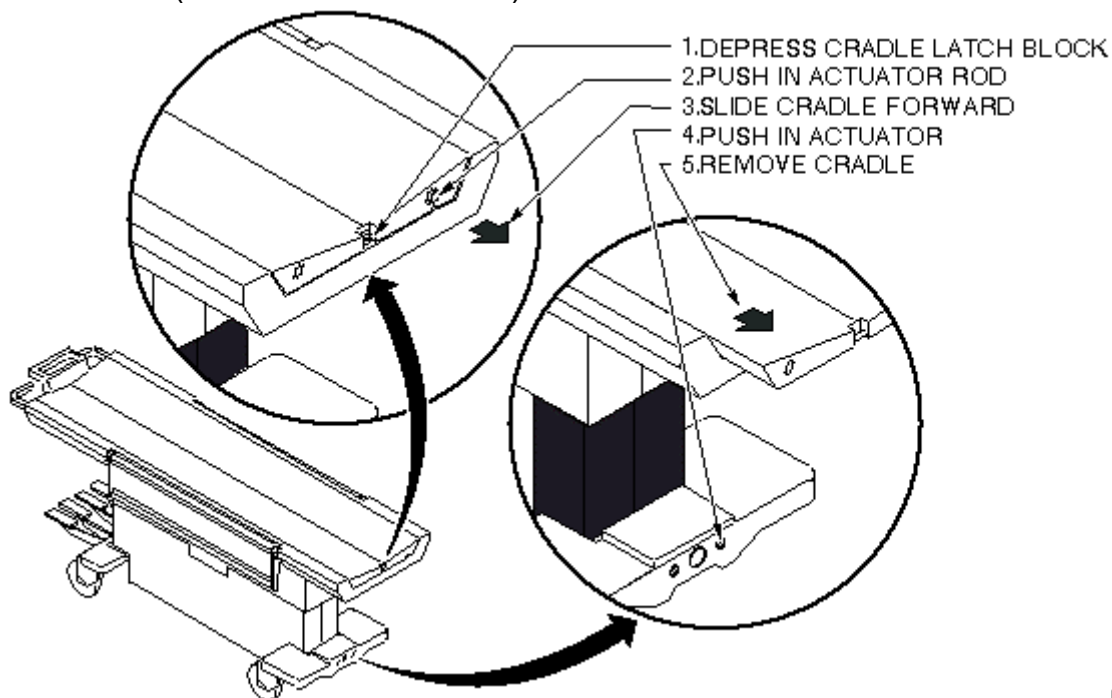

L2655A

PATIENT TRANSPORT CRADLE REMOVAL
ILLUSTRATION L2655A

5. Loosen set screw at interlock lever.
6. Unfasten cable end from interlock lever.
7. Pull cable out of pin stop on flipper bar (see Illustration L2675A).
8. Insert new cable through head of pin stop so that swaged end fits into hole in pin stop (see Illustration L2675A).

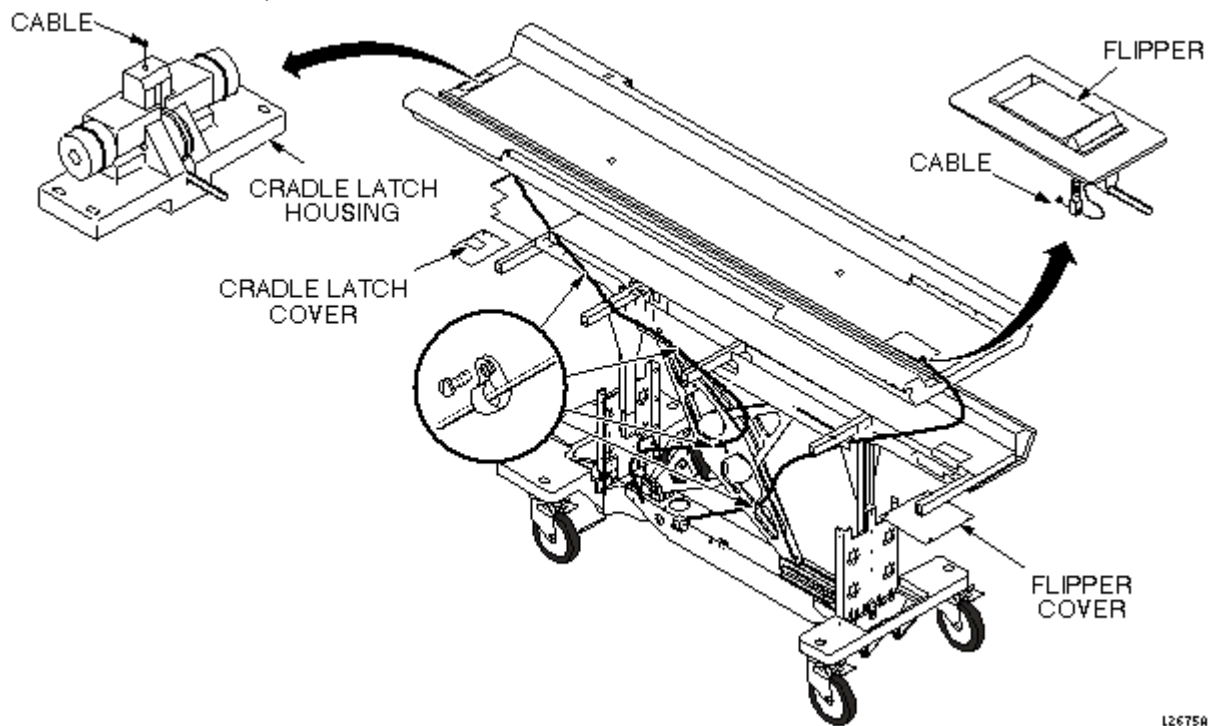


L2675A

INTERLOCK CABLE AND MAIN CRADLE LATCH CABLE ROUTING
ILLUSTRATION L2675A


9. Route cable through cable casing and insert through hole in interlock lever.

10. Loop cable twice through hole in interlock lever and secure with set screw. Tighten securely to prevent cable from slipping (see Illustration L2676A).

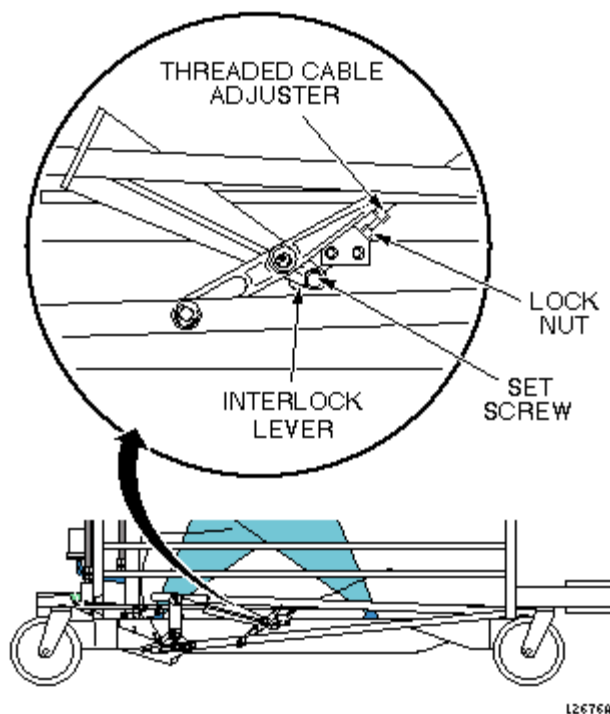

INTERLOCK LEVER
ILLUSTRATION L2676A

11. Adjust for indicated gap in down valve linkage when cradle is not in position. Use threaded cable adjuster for fine adjustments. Tighten lock nut after adjustment is completed, then cut off excess length of cable (see Illustration L2677A).

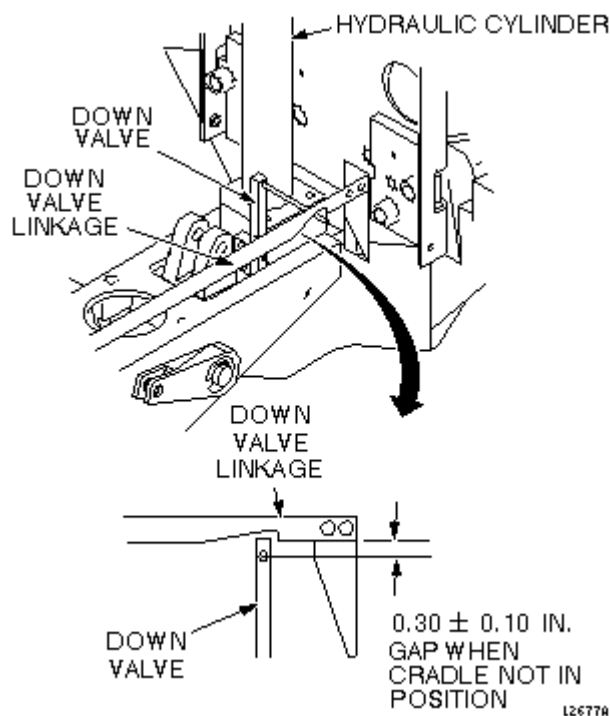


DOWN VALVE LINKAGE ADJUSTMENT
ILLUSTRATION L2677A

12. Replace cradle and upper and lower side covers. Note position of Signa logo on outer covers and position toward front of table (see Illustration L2654A).


11- MAIN CRADLE LATCH CABLE

1. Undock patient transport, and move it out of exam room.
2. Remove upper and lower side covers (see Illustration L2654A).

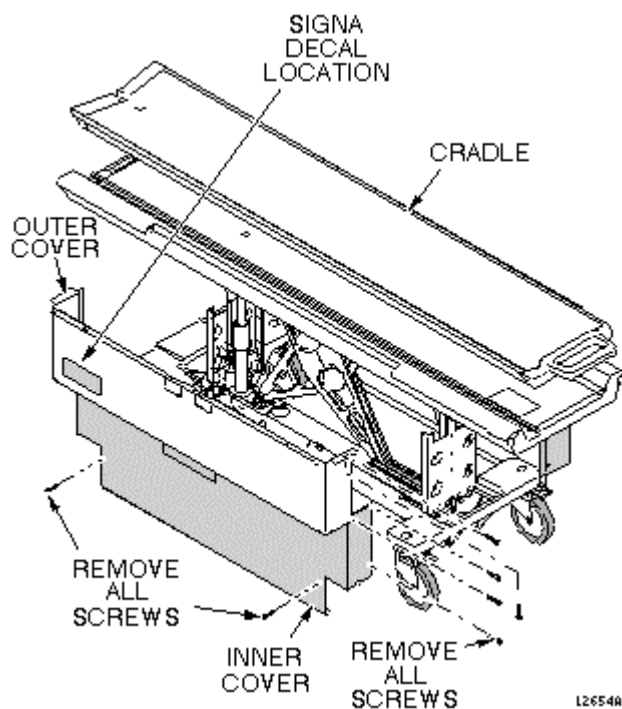

REMOVAL OF PATIENT TABLE COVERS
ILLUSTRATION L2654A


3. Remove cradle (see Illustration L2655A).

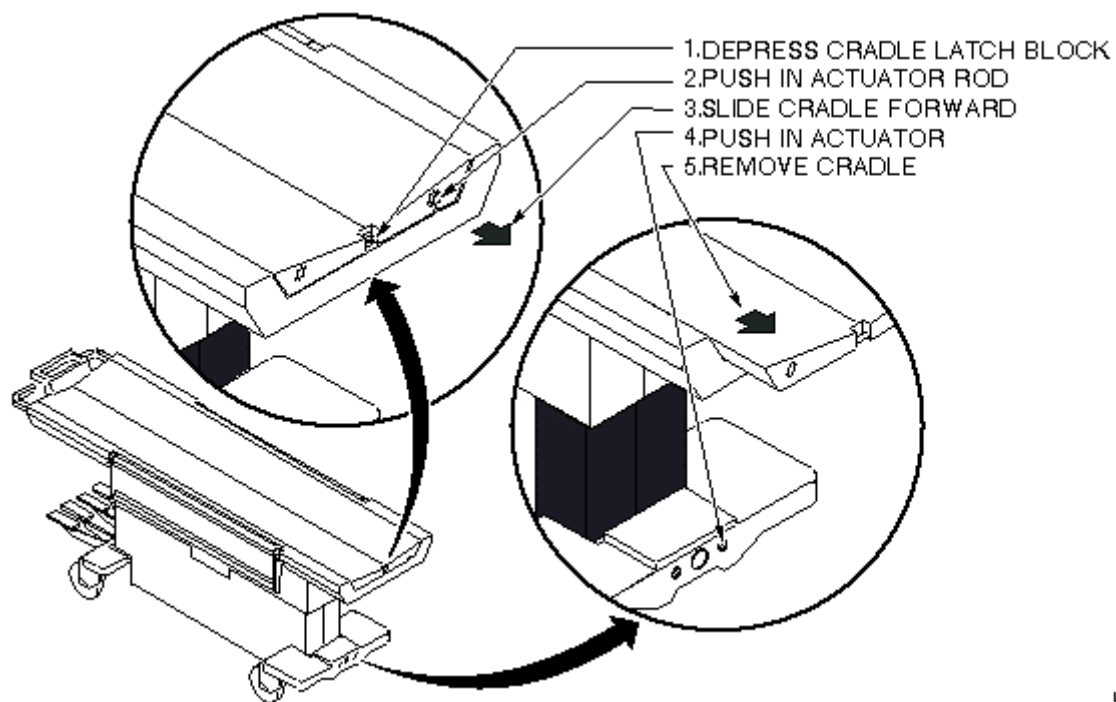

L2655A

PATIENT TRANSPORT CRADLE REMOVAL
ILLUSTRATION L2655A


4. Loosen set screw at down link clevis.
See Illustration L2678A).

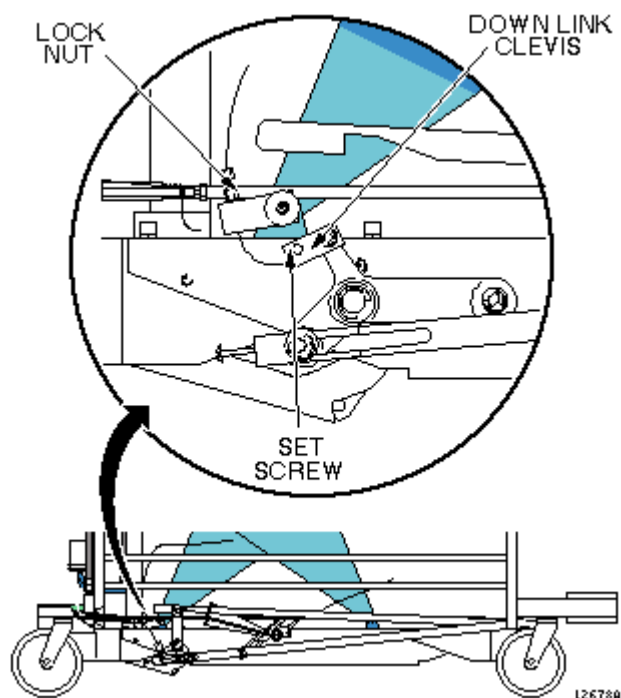

MAIN CRADLE LATCH DOWN LINK LEVER

ILLUSTRATION L2678A

5. Remove four bolts securing cradle latch housing and cradle latch cover (see Illustration L2675A).

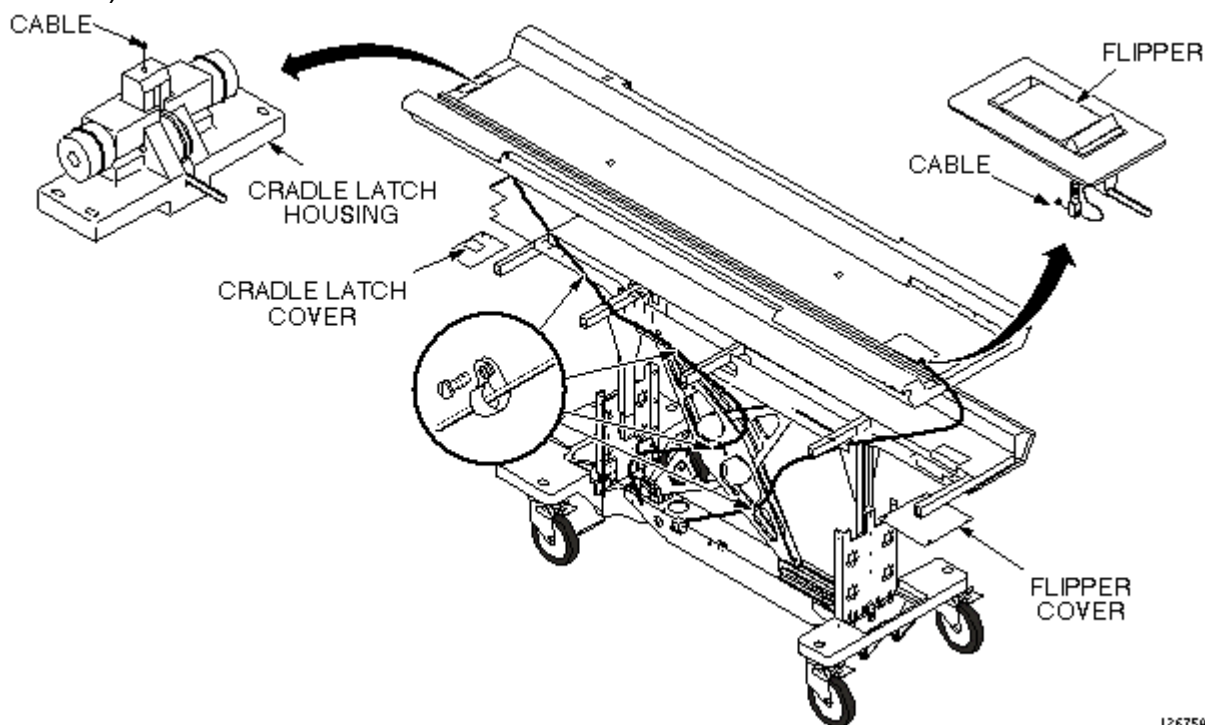


L2675A

INTERLOCK CABLE AND MAIN CRADLE LATCH CABLE ROUTING

ILLUSTRATION L2675A


6. Unfasten cable end from down link clevis (see Illustration L2678A).

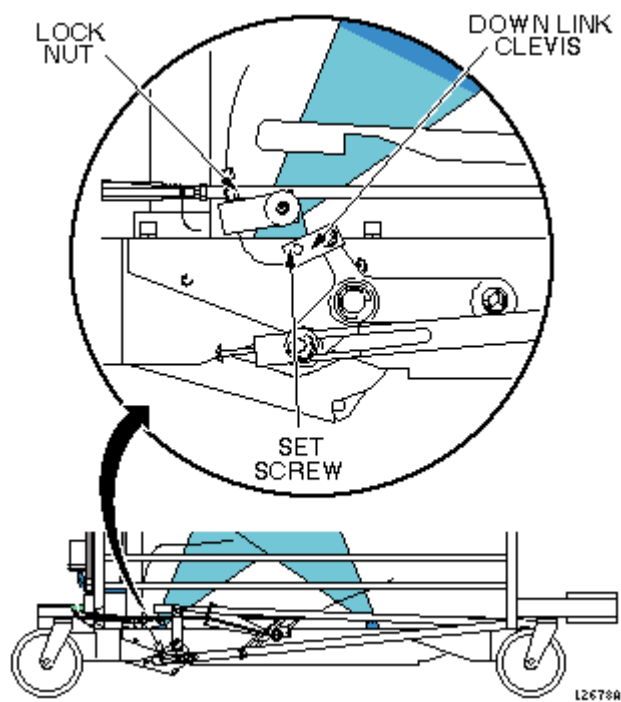

MAIN CRADLE LATCH DOWN LINK LEVER

ILLUSTRATION L2678A

7. Pull cable out of latch pin on front of patient transport (see Illustration L2675A).

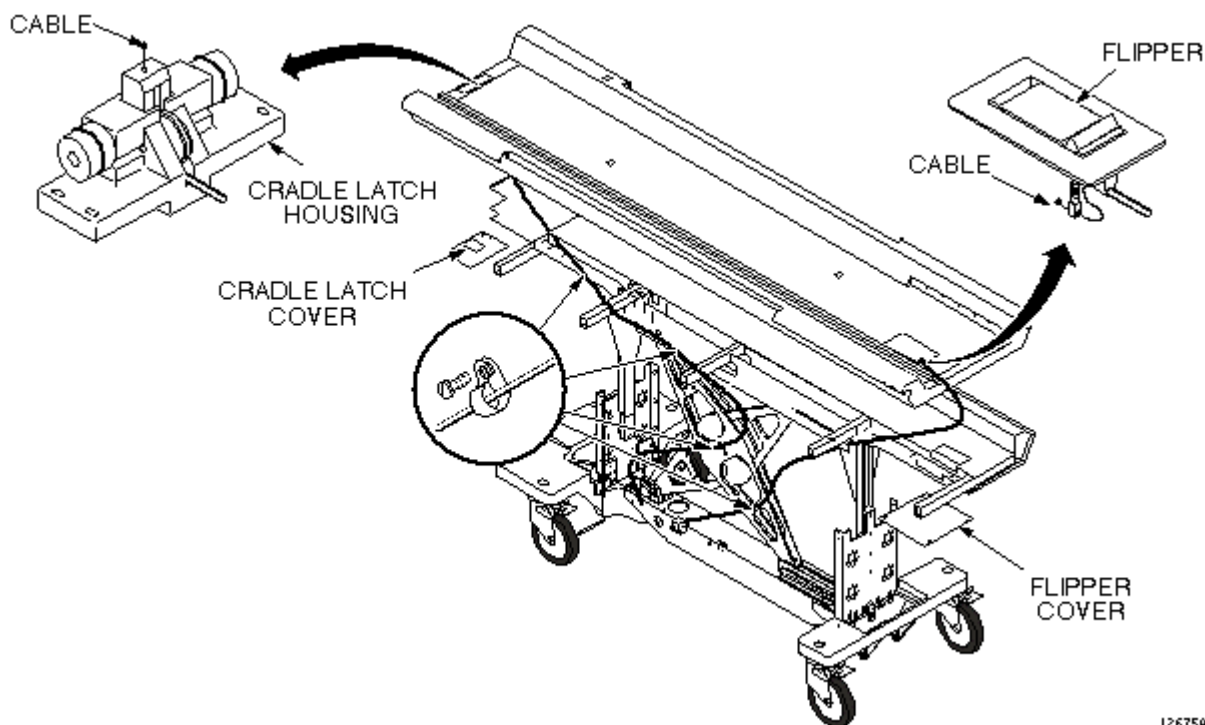


L2675A

INTERLOCK CABLE AND MAIN CRADLE LATCH CABLE ROUTING

ILLUSTRATION L2675A


8. Insert new cable through head of latch pin so that swaged end fits into hole in latch pin. See Illustration L2679A for details on installing cable.

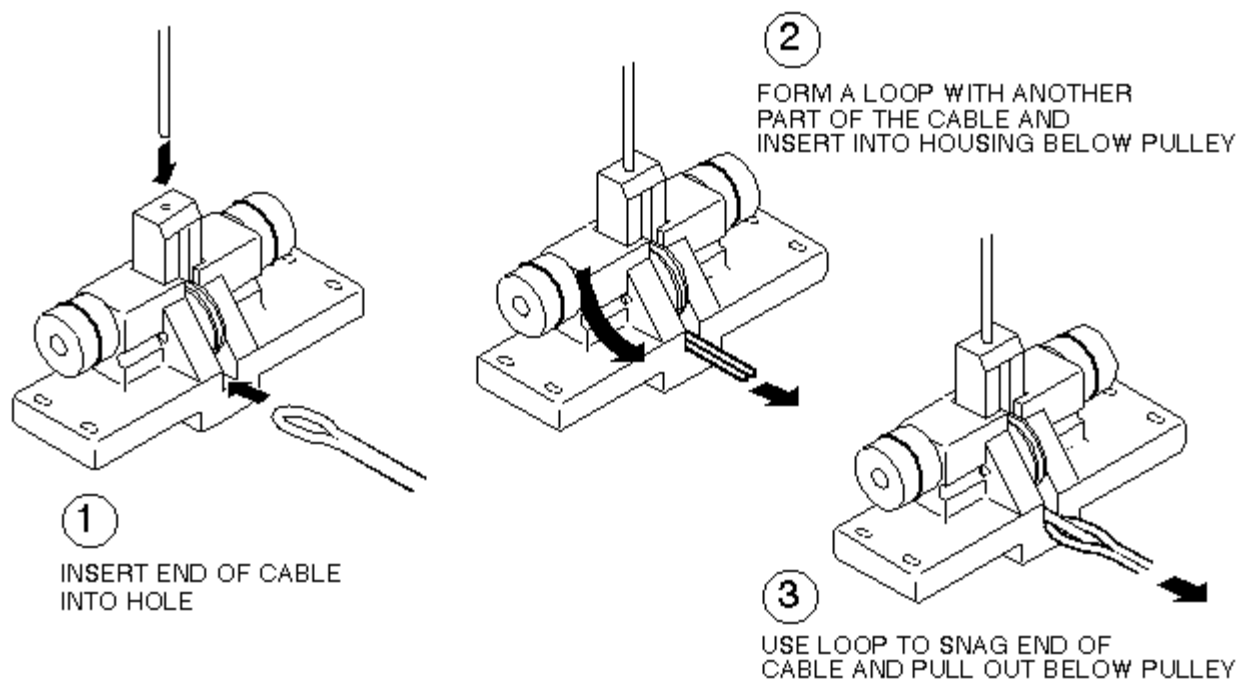

L2679A


INSTALLATION OF MAIN CRADLE LATCH CABLE
ILLUSTRATION L2679A

9. Route cable through cable casing and insert through cable clevis (see Illustration L2678A).
10. Loop cable twice through hole in down link clevis (see Illustration L2678A).
11. Tighten set screw (see Illustration L2678A).

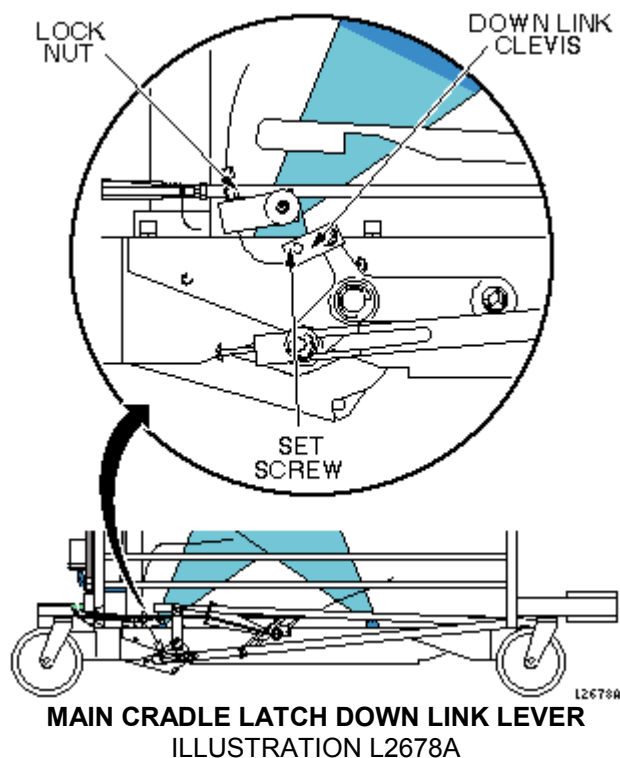


12. Return patient transport to exam room.
13. Dock patient transport.

14. Adjust cable so that when docked, main latch will be flush with casting. Refer to procedure for Lightweight Patient Transport, Section 6-2, Cradle Release Block Retraction Adjustment. Cut off excess length after adjustment is complete. Secure adjustment by tightening lock nut after adjustment is complete (see Illustration L2678A).

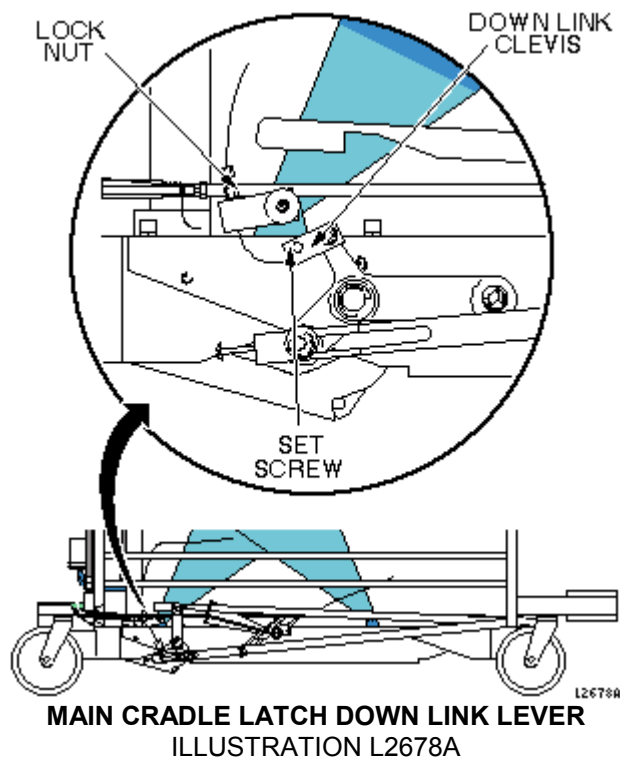

15. Replace cradle and upper and lower side covers. Note position of Signa logo on outer covers and position toward front of table (see Illustration L2654A).

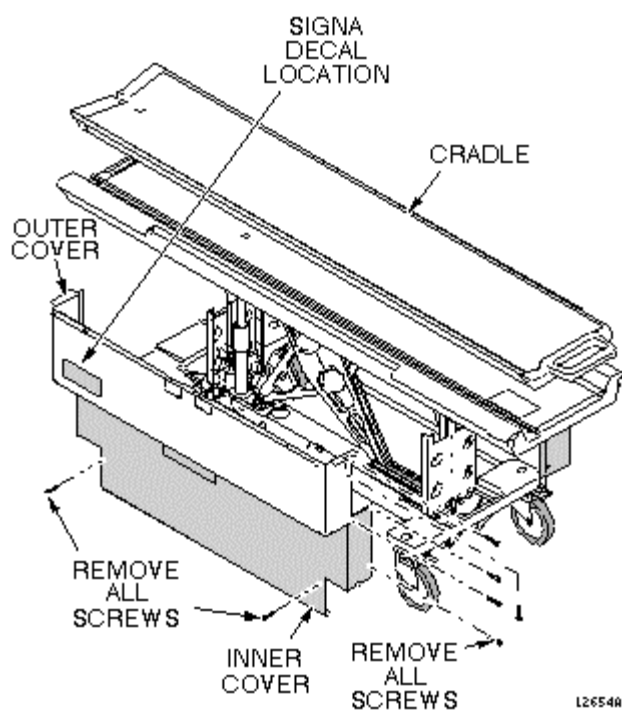

REMOVAL OF PATIENT TABLE COVERS
ILLUSTRATION L2654A

12- HYDRAULIC FILTER REPLACEMENT


Perform the following procedure once each year to maintain smooth, responsive operation of the patient transport hydraulics.

Required Materials:

Filter Kit	46-307878G1	qty. 1
Hydraulic Oil	46-230006P2	qty. 1
Absorbent Toweling		
13/16" open end wrench		qty 1
3/4" open end wrench		qty 2
Ear Syringe		qty 1


12-1 Remove Hydraulic Filter

1. Undock table and remove from magnet room.
2. Raise patient transport table to highest elevation with Up foot pedal.
3. Remove patient transport side covers on right side (when facing magnet enclosure).
4. Lower the safety latch by disconnecting spring from scissors assembly (see Illustration L2680A).

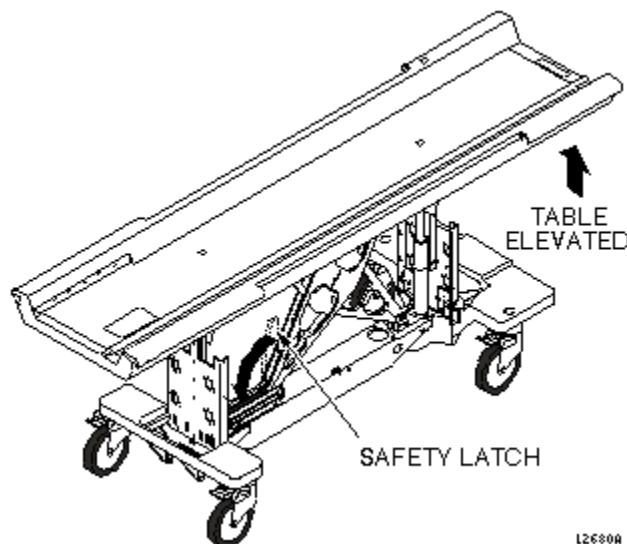

LOWERING SAFETY LATCH
ILLUSTRATION L2680A

5. Release hydraulic pressure by gently pressing Down foot pedal until table is supported by safety latch.
6. Place absorbent material or a drip pan on the floor under the cylinder block are to catch any dripping of oil.
7. Raise patient transport table with UP foot pedal (3 or 4 pumps) until table support bar is out of the way. This provides enough clearance to get a 13/16" wrench onto high pressure hose fitting (see Illustration L2681A).

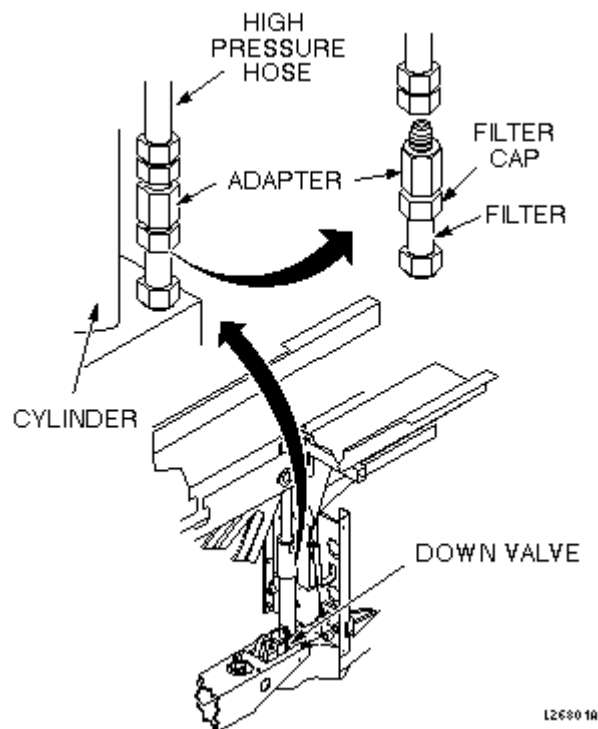


LOCATION OF HYDRAULIC FILTER
ILLUSTRATION L2681A

8. Turn wrench only enough to loosen the hose fitting initially, but not to let any oil leak out. Then lower the table using the Down foot pedal until table is again supported by safety latch.
9. Remove high-pressure hose from adapter. A small amount of oil will spill out of hose during removal; raise hose as high as possible to limit oil spillage.



Equipment damage possibility. Do not damage gasket with wrench during disassembly, if reusing old gasket. It is better to replace gasket with a new one, to prevent oil from leaking.

10. Remove filter cap from filter using two 3/4" wrenches (see Illustration L2682A).

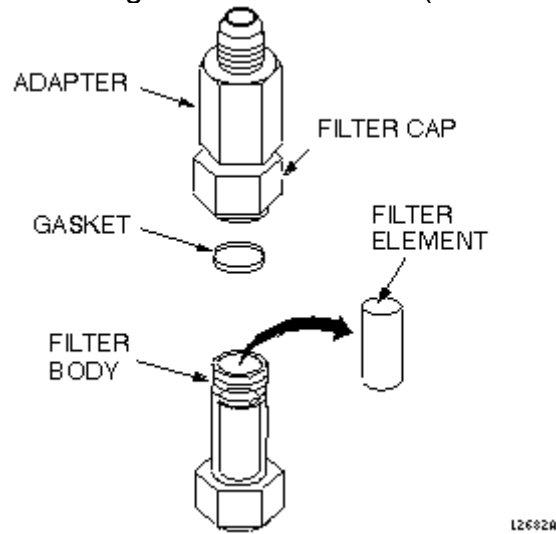

REMOVE FILTER ELEMENT
ILLUSTRATION L2682A

11. Remove dirty oil from top of filter.

12. Remove gasket and filter element from filter body and throw away.

12-2 Install Replacement Hydraulic Filter

1. Clean outside threads (both ends) of filter body with a clean cloth.

2. Clean interior threads of filter cap with a clean cloth.

3. Install new filter element
(46-230889P1), hole end down, in filter body.



Equipment damage possibility. When reassembling filter cap, use care not to damage gasket; ensure that gasket is properly seated before tightening with wrench. If gasket is deformed, oil will leak out of this connection when reassembled.


4. Install filter cap and new gasket
(46-230889P2) onto filter body. Ensure that new gasket seats properly before tightening with wrench.

5. Reconnect high-pressure hose to cylinder block. Tighten by hand, then raise the table with Up foot pedal to allow clearance for wrench; tighten hose fitting securely.

6. Raise table with several pumps on the Up foot pedal and check for oil leaks. If leaks occur, tighten connections snugly, but do not over tighten; brass fittings will deform if over tightened.
7. Proceed to step 5 of the next section.

13- BLEEDING PATIENT TRANSPORT HYDRAULICS


1. Unlock patient transport and remove from magnet room.
2. Raise patient transport table to highest elevation with Up foot pedal.
3. Remove patient transport side covers on right side (when facing magnet enclosure).
4. Lower the safety latch by disconnecting spring from scissors assembly (see Illustration L2680A).

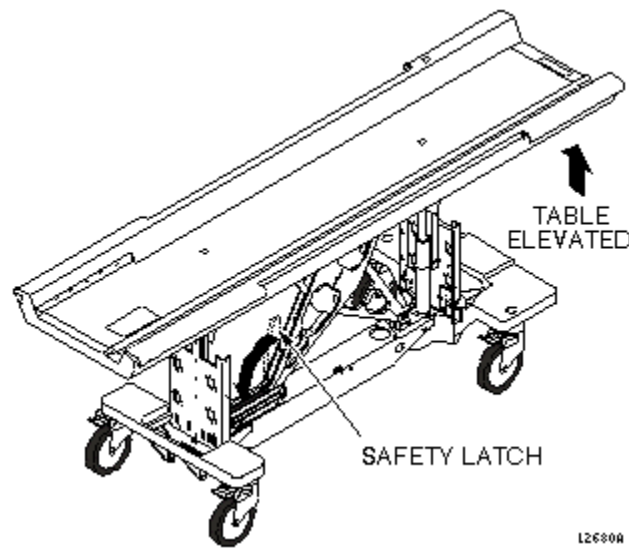

LOWERING SAFETY LATCH
ILLUSTRATION L2680A


5. Release hydraulic pressure by gently pressing Down foot pedal until table is supported by safety latch.

6. Hold the Down valve on the lift cylinder open while pumping the Up foot pedal at the same time. This can be done by one person by propping the Down valve in the open position using a standard blade screwdriver between the Down valve and the cylinder block. See Illustration L2681A for location of Down valve.

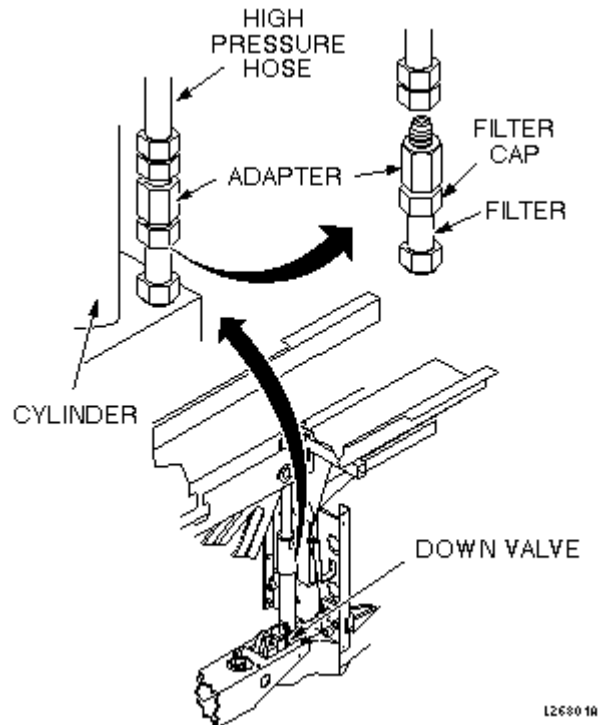

LOCATION OF HYDRAULIC FILTER
ILLUSTRATION L2681A

7. After all the air bubbles pass through the clear tubes (back to the reservoir), stop pumping and release the Down valve.
8. Remove the screwdriver to close the Down valve.
9. Purge the lift cylinder by pumping the table up to the top.
10. Open the bleed screw at the cylinder top a few turns until oil flows steadily from the screw after the air bubbles purge.
11. Tighten bleed screw (just snug).
12. With the table at highest elevation, raise the safety latch and reconnect spring to scissors assembly.
13. Lower table to its lowest position using Down foot pedal.
14. Fill table reservoir with hydraulic oil (46-230006P2) so that fluid is 1/4 to 1/2 inch above *full* line when table is at its lowest position (see Illustration L2593A).

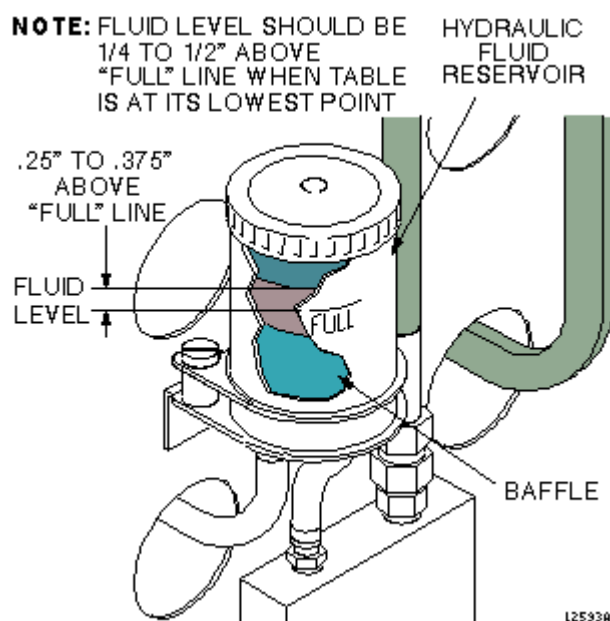

PATIENT TABLE RESERVOIR FILL LEVEL
ILLUSTRATION L2693A

15. Replace side covers.

14- SIGNA TABLE ARM BOARD BUMPER FIELD REPAIR/REPLACEMENT PROCEDURE

1- Tools and Materials

TABLE 14-1
TOOLS AND MATERIALS

Narrow (1 inch) Putty Knife	Local Purchase
Rubber Protector such as "Armor All" (or equivalent)	Local Purchase
Sharp Knife or Utility Knife	Local Purchase
Towels or Rags	Local Purchase
Adhesive Remover such as "GOO GONE" (or equivalent)	Local Purchase
Isopropyl Alcohol 46-183000P164	\$ 0.87 ICV *
Neoprene Gloves (100 count box) 46-194427P400 (lg) P401 (xlg)	\$ 12.76 ICV *
Adhesive (one quart ...fixes many bumpers) 46-170079P1	\$ 5.63 ICV *
Brush (for application of adhesive) 46-194427P218	\$ 2.00 ICV *
Bumper for Arm Board (if lost) 46-271017P2	\$ 5.74 ICV *

* Prices are approximant and given in USA dollars.

Approximant time needed to effect repairs as listed in this procedure is 54 minutes.

The creation of this procedure is the result of work done by Everest Castillo, Thomas Hulscher, Rick Neckerman, Kevin O'Connell, Robert Ralko, members of FIELD SERVICE SIX SIGMA GREEN BELT TEAM.

2- Repair/Replacement Procedure

1. Remove table from magnetic field to a ventilated location.
2. Remove partially attached bumper(s) (if present).
3. Clean arm board bumper channel(s) completely with alcohol and/or adhesive remover, putty knife and towel/rag. See Illustration 14-1.

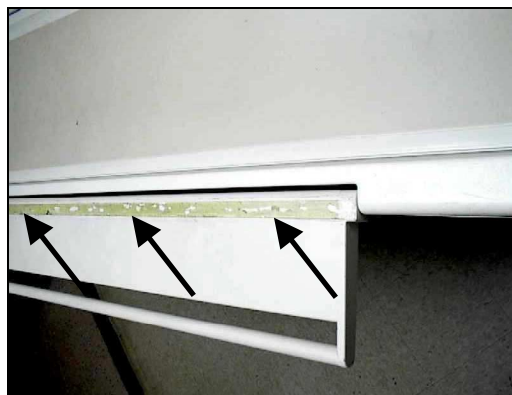


ILLUSTRATION 14-1
ARM BOARD BUMPER CHANNEL

4. Clean old glue off bumper(s) with alcohol or adhesive remover/clean. Or if installing new bumpers, clean new bumper(s) with alcohol.

5. Apply adhesive to bumper channel and bumper surface per manufacturer directions (use gloves and brush). Install bumper in channel with some over lap at each end. Apply pressure to entire bumper to insure good adhesion. Repeat on second bumper if necessary. See Illustration 14-2.

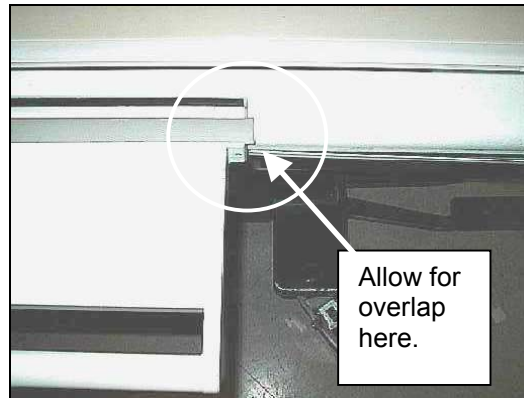


ILLUSTRATION 14-2
BUMPER INSTALLED WITH OVERLAP

6. Cut bumper ends at ~45 degree angles to match arm board edges. This will prevent catching on sheets, patients etc. See Illustration 14-3.



ILLUSTRATION 14-3
CUT BUMPER ENDS

7. Reduce surface friction by applying the rubber protector (Armor All if available) to bumper surface(s).
8. Return table to customer for use.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
A	April 17, 1998	W. Pasciak	Initial release converted to Word
0	April 18, 2001	Hawthorne	Added Section 14, Signa Table Arm Board Bumper Field Repair/Replacement Procedure. Also added the word "Repairs" to the title.