

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1 - OVERVIEW	2
2 - RESISTIVE SHIM SUPPLY REMOVAL	3
3 - RESISTIVE SHIM SUPPLY INSTALLATION	10
REVISION HISTORY	11

1 - OVERVIEW

This section applies to *TwinSpeed* Magnet systems and explains the procedure for the replacement of the Resistive Shim Supply.

Parts mentioned are listed in Table 1-1:

TABLE 1-1
FRU(1) PART LIST

Manufacturer's Part Number	GEMS Part Number	Item Description
	2257702	Resistive Shim Supply

Special tools needed to replace the Resistive Shim Supply are listed in Table 1-2:

TABLE 1-2
SPECIAL TOOLS NEEDED

Item	Description	GE Part Number
1	Universal Lift Hoist Kit	46-317266G5
2	Small hand tools including a Phillips-head screwdriver	
3	Two medium-sized crescent wrenches	

2 - RESISTIVE SHIM SUPPLY REMOVAL

DANGER!!


POSSIBLE PERSONAL INJURY! AVOID SERIOUS INJURY OR DEATH BY ELECTROCUTION. REMOVE POWER FROM THE TAC CABINET BEFORE ATTEMPTING TO REMOVE THE RESISTIVE SHIM SUPPLY. VERIFY THAT THE RFI/AMP AND SSM BREAKERS ON THE FRONT OF THE PDU ARE OFF. LOCK AND TAG OUT THE PLEXI-GLASS COVER ON THE FRONT OF THE PDU.

1. Lock out and tag out the TAC cabinet. (Refer to 8.X Service Methods CD-ROM, Safety (SYSMSA2.DOC), Section 6, OSHA LOCKOUT/TAGOUT REQUIREMENTS.)

DANGER!!

CABINET TIP HAZARD! POSSIBLE PERSONAL INJURY! CONFIRM THAT THE ALUMINUM ANTI-TIP LEGS STORED IN THE REAR OF THE CABINET ARE INSTALLED ON EITHER SIDE OF THE BASE OF THE CABINET BEFORE ATTEMPTING TO REMOVE AN RF AMPLIFIER.

2. Using the screws included with the anti-tip legs, secure the aluminum anti-tip legs to the threaded holes on each side of the base of the TAC Cabinet. See Illustration 2-1.

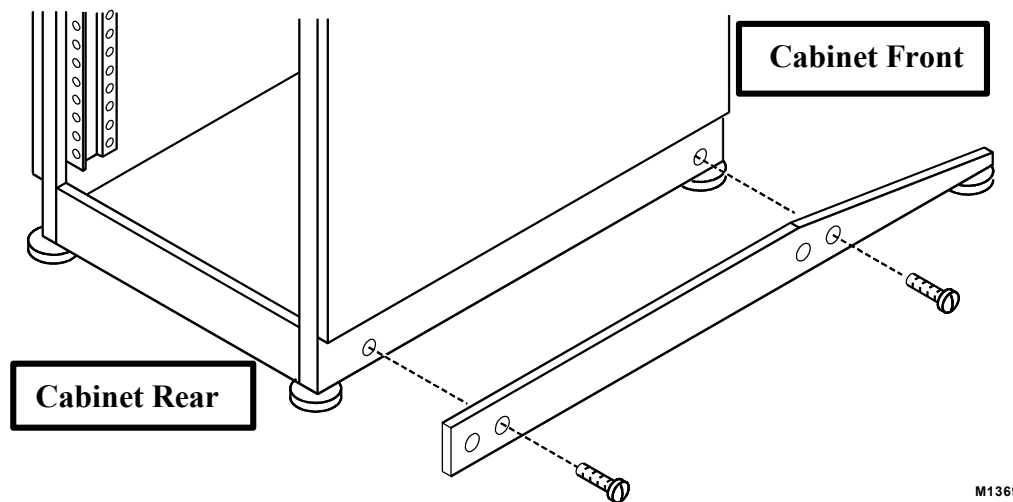

INSTALLING ANTI-TIP LEGS ON SIDES OF CABINET
ILLUSTRATION 2-1


3. Remove the front cover and the top bracket from the TAC cabinet.
4. Open the rear TAC cabinet door.
5. Remove the four (4) screws that attach the fan module to the TAC cabinet. Set the fan module in the cabinet.
6. Disconnect the data, COM 1, shim output, and power cables from the back of the resistive shim supply.

7. Locate and open the Universal Lift Hoist Kit.

NOTE


The Universal Lift Hoist is assembled from parts stored in the Universal Lift Hoist Kit.

8. Locate the three (3) large Steel Rail sections and place them on a flat surface for assembly. See Illustration 2-2.

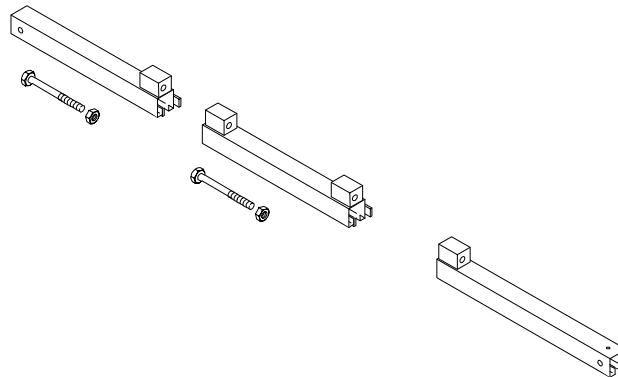

STEEL RAIL SECTIONS, HEX BOLTS, AND HEX NUTS
ILLUSTRATION 2-2

9. Arrange the rail sections in the proper order (rear, center, then front) and fasten them together with 4.75" X 1/2" caps crews and hex nuts.
10. Use the crescent wrenches to tighten the cap screws.
11. Assemble and fasten the elevating brackets to the rail assembly. See Illustration 2-3.

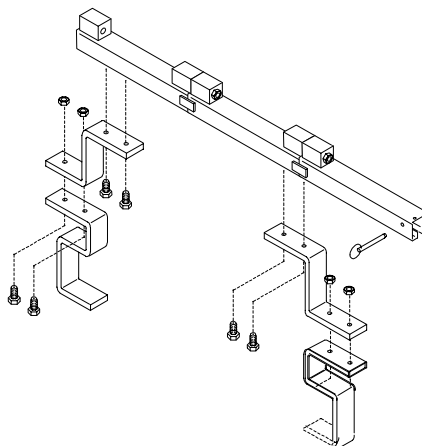

HOIST TOP RAIL PARTS WITH ELEVATING BRACKETS
ILLUSTRATION 2-3

12. Center the Steel Rail and elevating bracket assembly at the top of TAC Cabinet, so the hoist assembly will face the front of the cabinet. See Illustration 2-4.




HOIST INSTALLED ON TAC CABINET
ILLUSTRATION 2-4

13. Tighten the bolts on the Steel Rail to secure it to the top of the cabinet.



VERIFY THAT ALL HARDWARE IS SECURELY FASTENED; ALL BOLTS AND NUTS MUST BE TIGHT PRIOR TO ATTACHING THE HOIST ASSEMBLY TO THE LIFT KIT.

14. Locate the Hoist Assembly in the Lift Kit. See Illustration 2-5.

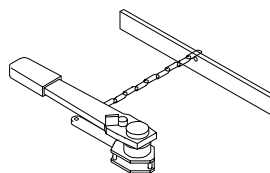

HOIST ASSEMBLY
ILLUSTRATION 2-5

15. Insert the Hoist Assembly into the open end of the front rail. Insert a locking pin into the holes near the front edge of the front rail. This locking pin prevents the Hoist Assembly from sliding out of the rails. See Illustration 2-6.




PLACEMENT OF HOIST ASSEMBLY INTO FRONT END OF STEEL RAIL ASSEMBLY
ILLUSTRATION 2-6

16. Remove the four (4) screws from the resistive shim supply.
17. Pull the resistive shim supply out of the cabinet until you see the STOP line. See Illustration 2-7.

DANGER!!

DO NOT PULL THE RESISTIVE SHIM SUPPLY OUT FARTHER THAN THE STOP LINE LOCATED ON THE TOP OF THE RESISTIVE SHIM SUPPLY.

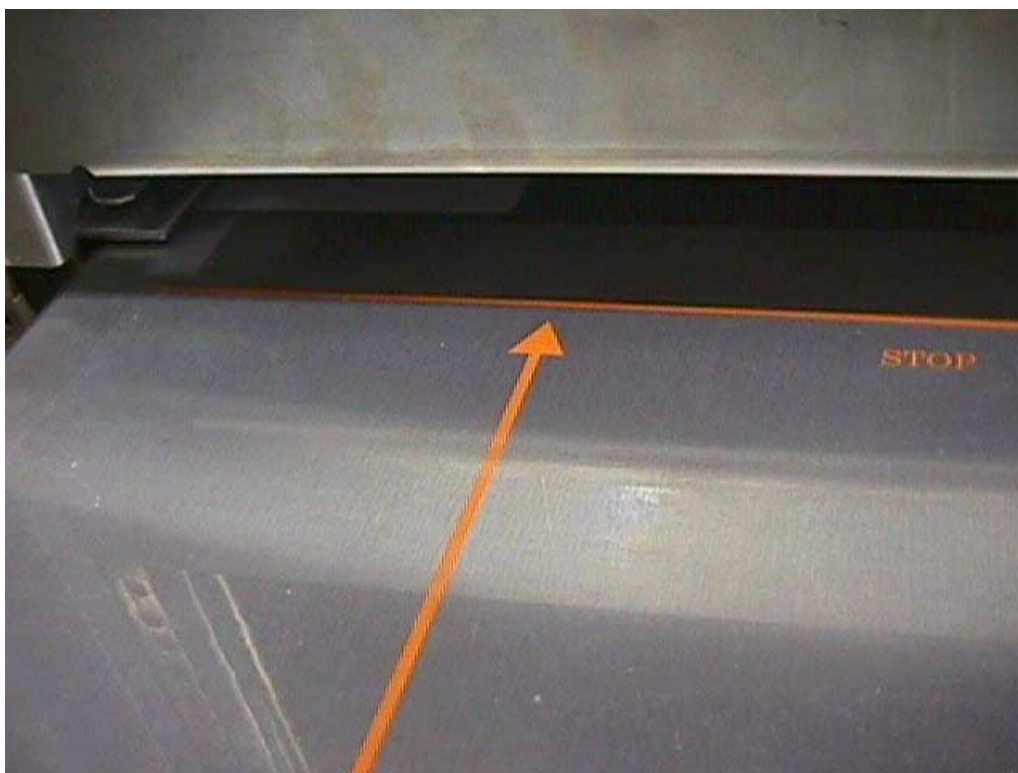

STOP LINE ON OF RESISTIVE SHIM SUPPLY
ILLUSTRATION 2-7

18. Attach the two side brackets to the resistive shim supply using the bolts on the side of the resistive shim supply. See Illustration 2-8.



HOIST SIDE BRACKETS
ILLUSTRATION 2-8

19. Attach the hoist bracket to the side brackets.
20. Hoist the resistive shim supply out of the cabinet and place it safely on the floor. See Illustration 2-9.



HOISTED RESISTIVE SHIM SUPPLY
ILLUSTRATION 2-9

21. Disconnect the hoist bracket from the side brackets and remove the side brackets from the resistive shim supply.

3 - RESISTIVE SHIM SUPPLY INSTALLATION

1. Attach the hoist side brackets to the resistive shim supply. See Illustration 2-8.
2. Attach the hoist bracket to the side brackets.
3. Hoist the resistive shim supply up and slide it part way into the TAC Cabinet.
4. Remove the hoist bracket and side brackets from the resistive shim supply.
5. Slide the resistive shim supply all the way into the cabinet.
6. Install the four (4) bolts that attach the resistive shim supply to the TAC Cabinet.
7. Attach the data, COM1, shim output, and power cables to the back of the resistive shim supply.
8. Remove the hoist assembly from front rail. See Illustration 2-6.
9. Loosen the bolts that secure the Steel Rail to the top of the cabinet.
10. Remove the Steel Rail from the cabinet.
11. Replace the top bracket and front cover for the TAC Cabinet.
12. Replace the fan module and secure with four (4) screws.
13. Close the rear TAC cabinet door.
14. Remove the anti-tip legs. See Illustration 2-1.
15. Remove the lock out and tag out for the TAC Cabinet.
16. Un-assemble and return the Universal Lift Hoist Kit components to the case.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	Oct. 24, 2001	K. Schraufnagel	Preliminary release