

GE Healthcare MR Class A eService Notes are copyright (c) protected to General Electric Company with all rights reserved.
Duplication and distribution of any class of MR eService Notes in any form, by entities other than GE Healthcare, is strictly prohibited.

SERVICE NOTE

11-JUN-2012

Class A

SN DOC1158554, Rev 1

1.5T HNS Latch FRU

EFFECTIVITY:

GE Healthcare Coils 1.5T HNS Coil M3335LG (2423192), M3335LH (2423192) and M7000AN (2423192).

PROBLEM:

In the past when the 1.5T HNS Coil had a latch failure it was required to return the complete coil. Returning the coil can be costly to the Service Organization.

SOLUTION:

Engineering released a Latch FRU kit (5365862). The new Latch FRU will repair the coil faster, more economically and prevent longer downtimes waiting for a replacement coil. The FRU kit contains one each of the anterior and posterior latch and screws if needed for replacement of the latch. The procedure for replacing the kit is contained within this document.

GE Healthcare MR Class A eService Notes are copyright (c) protected to General Electric Company with all rights reserved. Duplication and distribution of any class of MR eService Notes in any form, by entities other than GE Healthcare, is strictly prohibited.

1.5T HNS Latch Replacement

1. Overview

This process is used to replace the 1.5T HNS Latch. This procedure is limited to P/N 2416329 (FRU 2423192).

2. Preliminary Requirements

2.1. Personnel Requirements

Required Persons	Procedure
1	25 mins

2.2. Tools and Test Equipment

Item	Qty
Standard Tool Set – 3/16x4 Screw Driver	1

2.3. Replacement Part

Item	P/N	Qty
1.5T HNS Latch	5365862	1

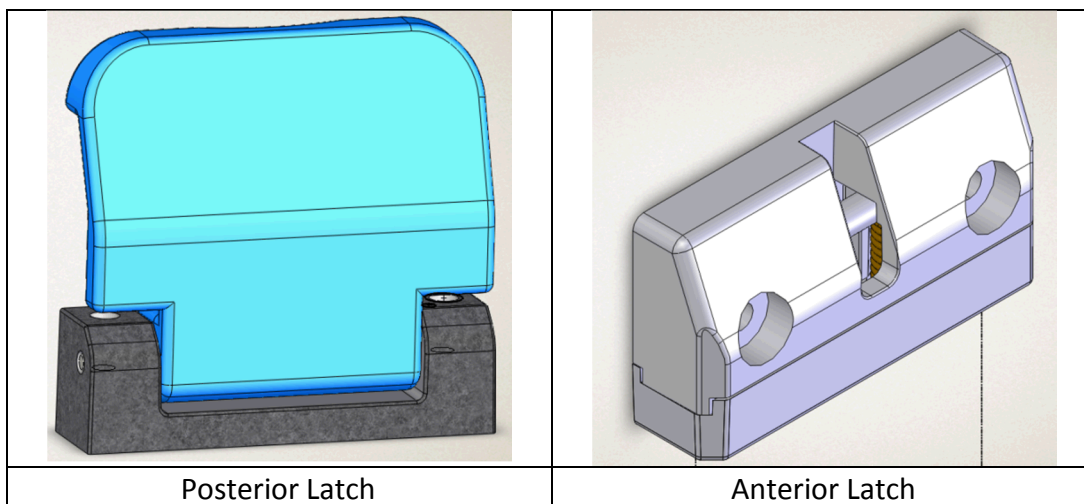


Figure 2.3

GE Healthcare MR Class A eService Notes are copyright (c) protected to General Electric Company with all rights reserved. Duplication and distribution of any class of MR eService Notes in any form, by entities other than GE Healthcare, is strictly prohibited.

3. Procedure

Note: Observe Anti-Static Precautions. Use a wrist strap and ESD workstation (ESD mat connected to ground and place to connect ESD wrist strap).

Note: Prior to completing the posterior latch repair it is recommended to complete an MCQA test.

Note: Repairs to the coil should take place outside the scan room.

- 3.1. If the posterior latch does not need replaced proceed to section 3.7
- 3.2. Position the 1.5T HNS Posterior Coil patient side down, exposing the screws. Remove the 16 screws as identified in Figure 3.1.

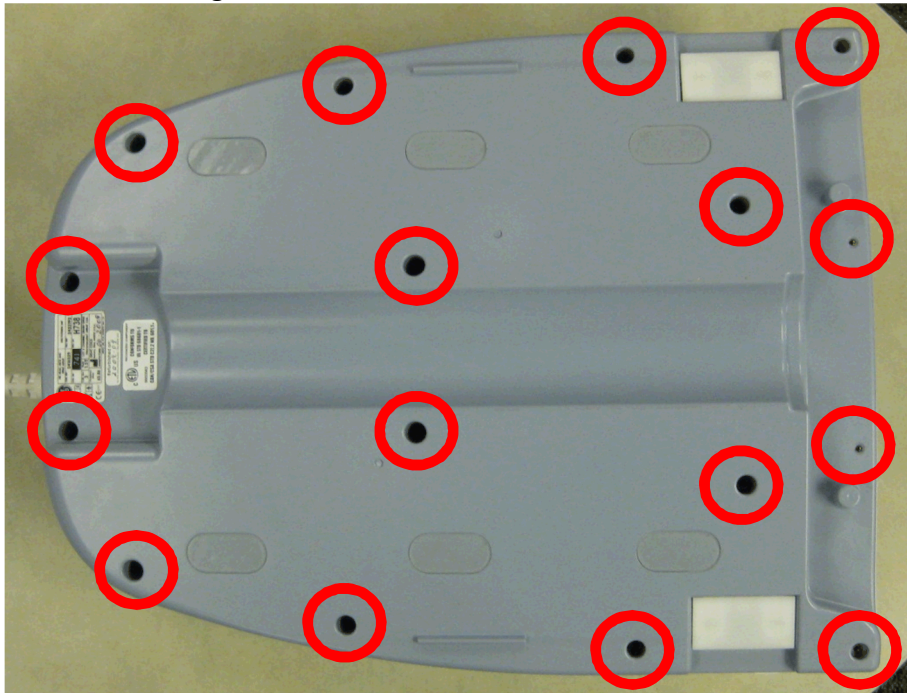


Figure 3.1

- 3.3. Remove the cover exposing the electronics of the coil as identified in Figure 3.2.

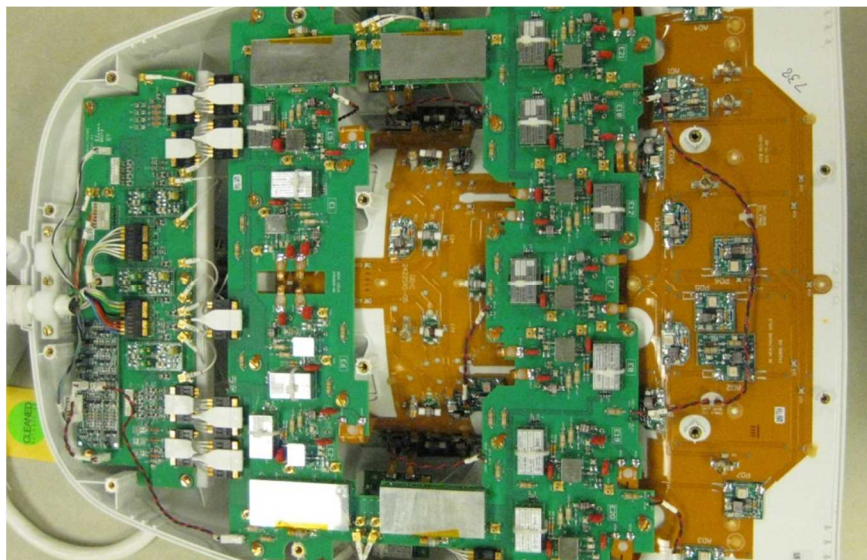


Figure 3.2

GE Healthcare MR Class A eService Notes are copyright (c) protected to General Electric Company with all rights reserved. Duplication and distribution of any class of MR eService Notes in any form, by entities other than GE Healthcare, is strictly prohibited.

- 3.4. Locate and remove the 2 screws connecting the posterior latch assembly and the former as seen in Figure 3.3 using the 3/16x4 Screw Driver.

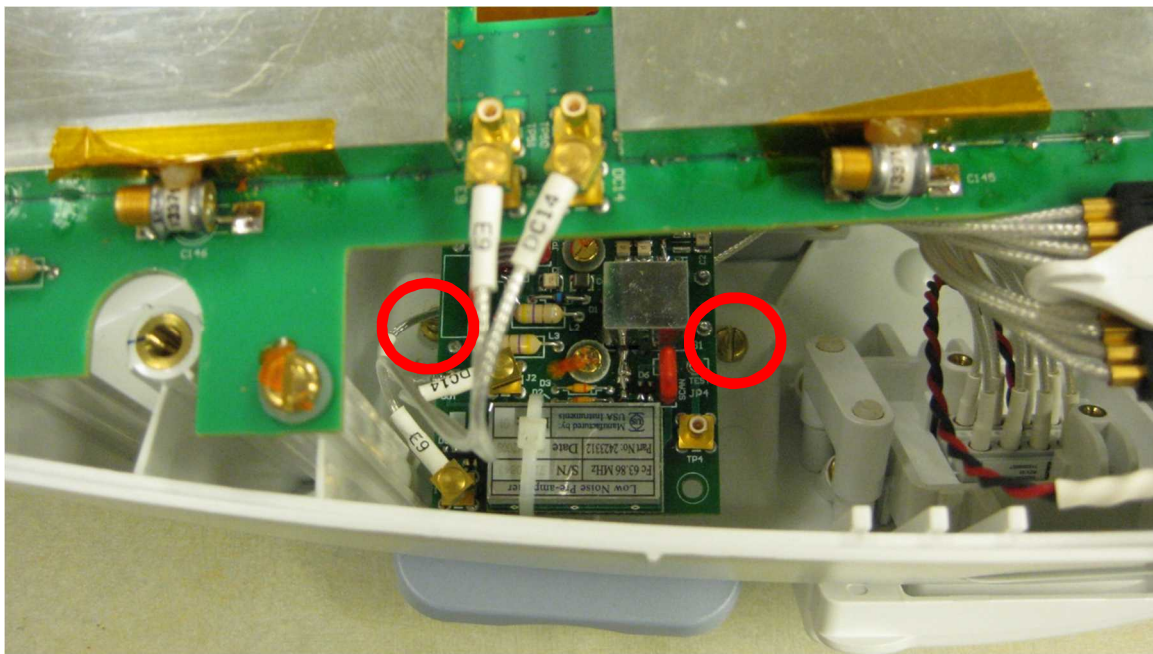


Figure 3.3

- 3.5. Install the new posterior latch assembly onto the posterior coil assembly with the 2 screws. Ensure the position of the latch is in the correct orientation prior to installing the screws as seen in Figure 3.4.

*Note: (Screws are part of the FRU kit for convenience. It is not necessary to replace them unless the screws are stripped or damaged.)

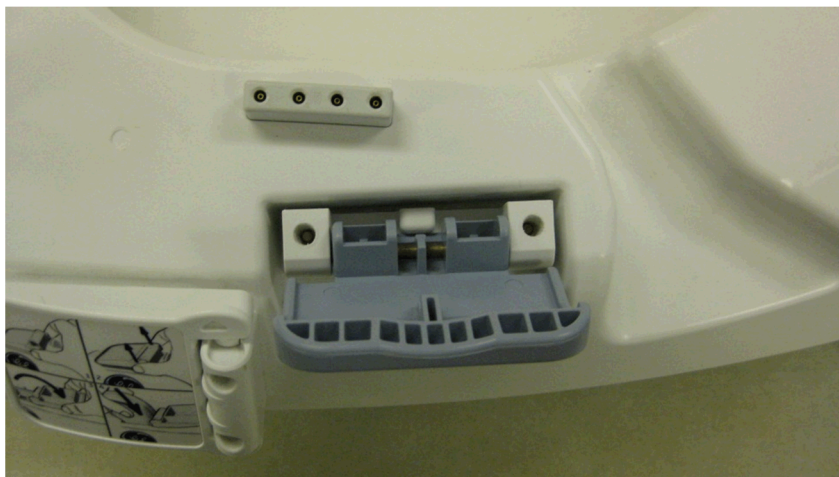


Figure 3.4

GE Healthcare MR Class A eService Notes are copyright (c) protected to General Electric Company with all rights reserved. Duplication and distribution of any class of MR eService Notes in any form, by entities other than GE Healthcare, is strictly prohibited.

- 3.6. Install the cover and 16 screws removed in step 3.1. If the anterior latch is not required to be replaced proceed to Section 4.
- 3.7. For replacing the anterior latch assembly, remove the two screws holding the latch clip to the 1.5T HNS anterior coil or the adapter block.
- 3.8. Place the anterior latch assembly in the correct orientation and attach it to the anterior coil or the adapter block with the 2 screws as seen in Figure 3.5.

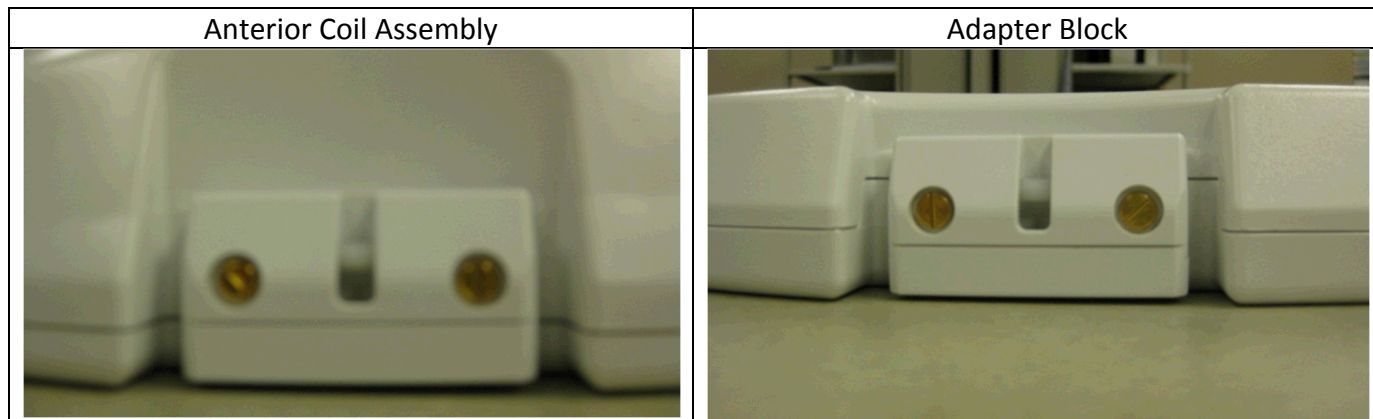


Figure 3.5

4. Finalization

- 4.1. Ensure the latch operates as intended without binding.
- 4.2. Perform a MCQA test to ensure the coil operates properly.

Author:

Stephen Scurria, Mechanical Engineer
GE Healthcare Coils
Aurora, OH

Contact:

Keith Richmond, LSI
GE Healthcare Coils
Aurora, OH
Ph: 330-995-5714