

GE Healthcare


MR750/MR450w OR-Compatible Patient Transport Table Service Manual



OPERATING DOCUMENTATION

5394275
Revision 1.0

WARNING

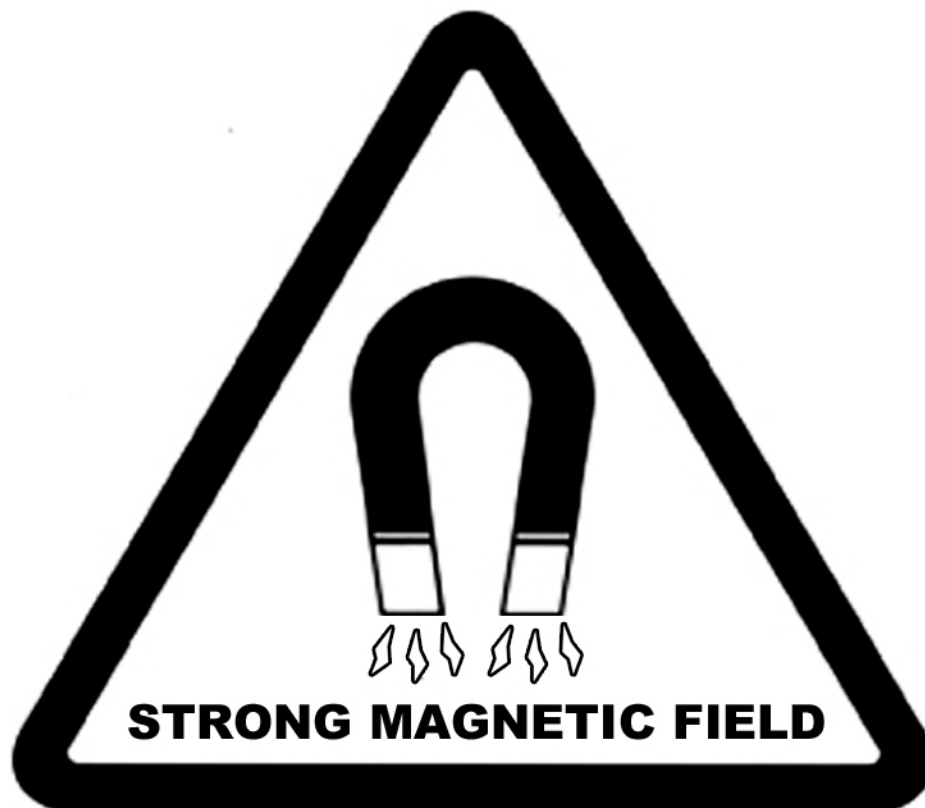

NO PACEMAKERS NO METALLIC IMPLANTS

Persons with pacemakers, neurostimulators or metallic implants must not enter the magnet area. Serious injury may result.



NO LOOSE METAL OBJECTS

Iron and steel materials must not be taken into the magnet area. Serious injury or property damage may result.

Important Information

LANGUAGE

ПРЕДУПРЕЖДЕНИЕ (BG)

Това упътване за работа е налично само на английски език.

- Ако доставчикът на услугата на клиента изиска друг език, задължение на клиента е да осигури превод.
- Не използвайте оборудването, преди да сте се консултирали и разбрали упътването за работа.
- Неспазването на това предупреждение може да доведе до нараняване на доставчика на услугата, оператора или пациента в резултат на токов удар, механична или друга опасност.

警告 (ZH-CN)

本维修手册仅提供英文版本。

- 如果客户的维修服务人员需要非英文版本，则客户需自行提供翻译服务。
- 未详细阅读和完全理解本维修手册之前，不得进行维修。
- 忽略本警告可能对维修服务人员、操作人员或患者造成电击、机械伤害或其他形式的伤害。

警告 (ZH-HK)

本服務手冊僅提供英文版本。

- 倘若客戶的服務供應商需要英文以外之服務手冊，客戶有責任提供翻譯服務。
- 除非已參閱本服務手冊及明白其內容，否則切勿嘗試維修設備。
- 不遵從本警告或會令服務供應商、網絡供應商或病人受到觸電、機械性或其他危險。

警告 (ZH-TW)

本維修手冊僅有英文版。

- 若客戶的維修廠商需要英文版以外的語言，應由客戶自行提供翻譯服務。
- 請勿試圖維修本設備，除非 您已查閱並瞭解本維修手冊。
- 若未留意本警告，可能導致維修廠商、操作員或病患因觸電、機械或其他危險而受傷。

UPOZORENJE (HR)

Ovaj servisni priručnik dostupan je na engleskom jeziku.

- Ako davatelj usluge klijenta treba neki drugi jezik, klijent je dužan osigurati prijevod.
- Ne pokušavajte servisirati opremu ako niste u potpunosti pročitali i razumjeli ovaj servisni priručnik.
- Zanimarite li ovo upozorenje, može doći do ozljede davatelja usluge, operatera ili pacijenta uslijed strujnog udara, mehaničkih ili drugih rizika.

**VÝSTRAHA
(CS)**

Tento provozní návod existuje pouze v anglickém jazyce.

- V případě, že externí služba zákazníkům potřebuje návod v jiném jazyce, je zajištění překladu do odpovídajícího jazyka úkolem zákazníka.
- Nesnažte se o údržbu tohoto zařízení, aniž byste si přečetli tento provozní návod a pochopili jeho obsah.
- V případě nedodržování této výstrahy může dojít k poranění pracovníka prodejního servisu, obslužného personálu nebo pacientů vlivem elektrického proudu, respektive vlivem mechanických či jiných rizik.

**ADVARSEL
(DA)**

Denne servicemanual findes kun på engelsk.

- Hvis en kundes tekniker har brug for et andet sprog end engelsk, er det kundens ansvar at sørge for oversættelse.
- Forsøg ikke at servicere udstyret uden at læse og forstå denne servicemanual.
- Manglende overholdelse af denne advarsel kan medføre skade på grund af elektrisk stød, mekanisk eller anden fare for teknikeren, operatøren eller patienten.

**WAARSCHUWING
(NL)**

Deze onderhoudshandleiding is enkel in het Engels verkrijgbaar.

- Als het onderhoudspersoneel een andere taal vereist, dan is de klant verantwoordelijk voor de vertaling ervan.
- Probeer de apparatuur niet te onderhouden alvorens deze onderhoudshandleiding werd geraadpleegd en begrepen is.
- Indien deze waarschuwing niet wordt opgevolgd, zou het onderhoudspersoneel, de operator of een patiënt gewond kunnen raken als gevolg van een elektrische schok, mechanische of andere gevaren.

**WARNING
(EN)**

This service manual is available in English only.

- If a customer's service provider requires a language other than English, it is the customer's responsibility to provide translation services.
- Do not attempt to service the equipment unless this service manual has been consulted and is understood.
- Failure to heed this warning may result in injury to the service provider, operator or patient from electric shock, mechanical or other hazards.

**HOIATUS
(ET)**

See teenindusjuhend on saadaval ainult inglise keeles.

- Kui klienditeeninduse osutaja nõuab juhendit inglise keelest erinevas keeles, vastutab klient tõlketeenuse osutamise eest.
- Ärge üritage seadmeid teenindada enne eelnevalt käesoleva teenindusjuhendiga tutvumist ja sellest aru saamist.
- Käesoleva hoiatuse eiramine võib põhjustada teenuseosutaja, operaatori või patsiendi vigastamist elektrilöögi, mehaanilise või muu ohu tagajärjel.

**VAROITUS
(FI)**

Tämä huolto-ohje on saatavilla vain englanniksi.

- Jos asiakkaan huoltohenkilöstö vaatii muuta kuin englanninkielistä materiaalia, tarvittavan käännöksen hankkiminen on asiakkaan vastuulla.
- Älä yritä korjata laitteistoa ennen kuin olet varmasti lukenut ja ymmärtänyt tämän huolto-ohjeen.
- Mikäli tätä varoitusta ei noudateta, seurauksena voi olla huoltohenkilöstön, laitteiston käyttäjän tai potilaan vahingoittuminen sähköiskun, mekaanisen vian tai muun vaaratilanteen vuoksi.

**ATTENTION
(FR)**

Ce manuel d'installation et de maintenance est disponible uniquement en anglais.

- Si le technicien d'un client a besoin de ce manuel dans une langue autre que l'anglais, il incombe au client de le faire traduire.
- Ne pas tenter d'intervenir sur les équipements tant que ce manuel d'installation et de maintenance n'a pas été consulté et compris.
- Le non-respect de cet avertissement peut entraîner chez le technicien, l'opérateur ou le patient des blessures dues à des dangers électriques, mécaniques ou autres.

**WARNUNG
(DE)**

Diese Serviceanleitung existiert nur in englischer Sprache.

- Falls ein fremder Kundendienst eine andere Sprache benötigt, ist es Aufgabe des Kunden für eine entsprechende Übersetzung zu sorgen.
- Versuchen Sie nicht diese Anlage zu warten, ohne diese Serviceanleitung gelesen und verstanden zu haben.
- Wird diese Warnung nicht beachtet, so kann es zu Verletzungen des Kundendienst-technikers, des Bedieners oder des Patienten durch Stromschläge, mechanische oder sonstige Gefahren kommen.

**ΠΡΟΕΙΔΟΠΟΙΗΣΗ
(EL)**

Το παρόν εγχειρίδιο σέρβις διατίθεται μόνο στα αγγλικά.

- Εάν ο τεχνικός σέρβις ενός πελάτη απαιτεί το παρόν εγχειρίδιο σε γλώσσα εκτός των αγγλικών, αποτελεί ευθύνη του πελάτη να παρέχει τις υπηρεσίες μετάφρασης.
- Μην επιχειρήσετε την εκτέλεση εργασιών σέρβις στον εξοπλισμό αν δεν έχετε συμβουλευτεί και κατανοήσει το παρόν εγχειρίδιο σέρβις.
- Αν δεν προσέξετε την προειδοποίηση αυτή, ενδέχεται να προκληθεί τραυματισμός στον τεχνικό σέρβις, στο χειριστή ή στον ασθενή από ηλεκτροπληξία, μηχανικούς ή άλλους κινδύνους.

**FIGYELMEZTETÉS
(HU)**

Ezen karbantartási kézikönyv kizárólag angol nyelven érhető el.

- Ha a vevő szolgáltatója angoltól eltérő nyelvre tart igényt, akkor a vevő felelőssége a fordítás elkészíttetése.
- Ne próbálja elkezdni használni a berendezést, amíg a karbantartási kézikönyvben leírtakat nem értelmezték.
- Ezen figyelmeztetés figyelmen kívül hagyása a szolgáltató, működtető vagy a beteg áramütés, mechanikai vagy egyéb veszélyhelyzet miatti sérülését eredményezheti.

**ADVÖRUN
(IS)**

Þessi þjónustuhandbók er aðeins fánleg á ensku.

- Ef að þjónustuveitandi viðskiptamanns þarfnast annas tungumáls en ensku, er það skylda viðskiptamanns að skaffa tungumálaþjónustu.
- Reynið ekki að afgreiða tækið nema að þessi þjónustuhandbók hefur verið skoðuð og skilin.
- Brot á sinna þessari aðvörun getur leitt til meiðsla á þjónustuveitanda, stjórnanda eða sjúklings frá raflosti, vélrænu eða öðrum áhættum.

**AVVERTENZA
(IT)**

Il presente manuale di manutenzione è disponibile soltanto in lingua inglese.

- Se un addetto alla manutenzione richiede il manuale in una lingua diversa, il cliente è tenuto a provvedere direttamente alla traduzione.
- Procedere alla manutenzione dell'apparecchiatura solo dopo aver consultato il presente manuale ed averne compreso il contenuto.
- Il mancato rispetto della presente avvertenza potrebbe causare lesioni all'addetto alla manutenzione, all'operatore o ai pazienti provocate da scosse elettriche, urti meccanici o altri rischi.

**警告
(JA)**

このサービスマニュアルには英語版しかありません。

- サービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。
- このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないでください。
- この警告に従わない場合、サービスを担当される方、操作員あるいは患者さんが、感電や機械的又はその他の危険により負傷する可能性があります。

**경고
(KO)**

본 서비스 매뉴얼은 영어로만 이용하실 수 있습니다.

- 고객의 서비스 제공자가 영어 이외의 언어를 요구할 경우, 번역 서비스를 제공하는 것은 고객의 책임입니다.
- 본 서비스 매뉴얼을 참조하여 숙지하지 않은 이상 해당 장비를 수리하려고 시도하지 마십시오.
- 본 경고 사항에 유의하지 않으면 전기 쇼크, 기계적 위험, 또는 기타 위험으로 인해 서비스 제공자, 사용자 또는 환자에게 부상을 입힐 수 있습니다.

**BRĪDINĀJUMS
(LV)**

Šī apkopes rokasgrāmata ir pieejama tikai angļu valodā.

- Ja klienta apkopes sniedzējam nepieciešama informācija citā valodā, klienta pienākums ir nodrošināt tulkojumu.
- Neveiciet aprīkojuma apkopi bez apkopes rokasgrāmatas izlasīšanas un saprašanas.
- Šī brīdinājuma neievērošanas rezultātā var rasties elektriskās strāvas trieciena, mehānisku vai citu faktoru izraisītu traumų risks apkopes sniedzējam, operatoram vai pacientam.

**ĮSPĖJIMAS
(LT)**

Šis eksploatavimo vadovas yra tik anglų kalba.

- Jei kliento paslaugų tiekėjas reikalauja vadovo kita kalba – ne anglų, suteikti vertimo paslaugas privalo klientas.
- Nemėginkite atlikti įrangos techninės priežiūros, jei neperskaitėte ar nesupratote šio eksploatavimo vadovo.
- Jei nepaisysite šio įspėjimo, galimi paslaugų tiekėjo, operatoriaus ar paciento sužalojimai dėl elektros šoko, mechaninių ar kitų pavojų.

**ADVARSEL
(NO)**

Denne servicehåndboken finnes bare på engelsk.

- Hvis kundens serviceleverandør har bruk for et annet språk, er det kundens ansvar å sørge for oversettelse.
- Ikke forsøk å reparere utstyret uten at denne servicehåndboken er lest og forstått.
- Manglende hensyn til denne advarselen kan føre til at serviceleverandøren, operatøren eller pasienten skades på grunn av elektrisk støt, mekaniske eller andre farer.

**OSTRZEŻENIE
(PL)**

Niniejszy podręcznik serwisowy dostępny jest jedynie w języku angielskim.

- Jeśli serwisant klienta wymaga języka innego niż angielski, zapewnienie usługi tłumaczenia jest obowiązkiem klienta.
- Nie próbować serwisować urządzenia bez zapoznania się z niniejszym podręcznikiem serwisowym i zrozumienia go.
- Niezastosowanie się do tego ostrzeżenia może doprowadzić do obrażeń serwisanta, operatora lub pacjenta w wyniku porażenia prądem elektrycznym, zagrożenia mechanicznego bądź innego.

**ATENÇÃO
(PT-BR)**

Este manual de assistência técnica encontra-se disponível unicamente em inglês.

- Se outro serviço de assistência técnica solicitar a tradução deste manual, caberá ao cliente fornecer os serviços de tradução.
- Não tente reparar o equipamento sem ter consultado e compreendido este manual de assistência técnica.
- A não observância deste aviso pode ocasionar ferimentos no técnico, operador ou paciente decorrentes de choques elétricos, mecânicos ou outros.

**ATENÇÃO
(PT-PT)**

Este manual de assistência técnica só se encontra disponível em inglês.

- Se qualquer outro serviço de assistência técnica solicitar este manual noutro idioma, é da responsabilidade do cliente fornecer os serviços de tradução.
- Não tente reparar o equipamento sem ter consultado e compreendido este manual de assistência técnica.
- O não cumprimento deste aviso pode colocar em perigo a segurança do técnico, do operador ou do paciente devido a choques eléctricos, mecânicos ou outros.

**ATENȚIE
(RO)**

Acest manual de service este disponibil doar în limba engleză.

- Dacă un furnizor de servicii pentru clienți necesită o altă limbă decât cea engleză, este de datoria clientului să furnizeze o traducere.
- Nu încercați să reparați echipamentul decât ulterior consultării și înțelegerii acestui manual de service.
- Ignorarea acestui avertisment ar putea duce la rănirea depanatorului, operatorului sau pacientului în urma pericolelor de electrocutare, mecanice sau de altă natură.

**ОСТОРОЖНО!
(RU)**

Данное руководство по техническому обслуживанию представлено только на английском языке.

- Если сервисному персоналу клиента необходимо руководство не на английском, а на каком-то другом языке, клиенту следует самостоятельно обеспечить перевод.
- Перед техническим обслуживанием оборудования обязательно обратитесь к данному руководству и поймите изложенные в нем сведения.
- Несоблюдение требований данного предупреждения может привести к тому, что специалист по техобслуживанию, оператор или пациент получит удар электрическим током, механическую травму или другое повреждение.

**UPOZORENJE
(SR)**

Ovo servisno uputstvo je dostupno samo na engleskom jeziku.

- Ako klijentov serviser zahteva neki drugi jezik, klijent je dužan da obezbedi prevodi-lačke usluge.
- Ne pokušavajte da opravite uređaj ako niste pročitali i razumeli ovo servisno uputstvo.
- Zanemarivanje ovog upozorenja može dovesti do povređivanja serviser, rukovaoca ili pacijenta usled strujnog udara ili mehaničkih i drugih opasnosti.

**UPOZORNENIE
(SK)**

Tento návod na obsluhu je k dispozícii len v angličtine.

- Ak zákazník poskytovateľ služieb vyžaduje iný jazyk ako angličtinu, poskytnutie prekladateľských služieb je zodpovednosťou zákazníka.
- Nepokúšajte sa o obsluhu zariadenia, kým si neprečítate návod na obsluhu a neporozumiete mu.
- Zanedbanie tohto upozornenia môže spôsobiť zranenie poskytovateľa služieb, obsluhujúcej osoby alebo pacienta elektrickým prúdom, mechanické alebo iné ohrozenie.

**ATENCION
(ES)**

Este manual de servicio sólo existe en inglés.

- Si el encargado de mantenimiento de un cliente necesita un idioma que no sea el inglés, el cliente deberá encargarse de la traducción del manual.
- No se deberá dar servicio técnico al equipo, sin haber consultado y comprendido este manual de servicio.
- La no observancia del presente aviso puede dar lugar a que el proveedor de servicios, el operador o el paciente sufran lesiones provocadas por causas eléctricas, mecánicas o de otra naturaleza.

**VARNING
(SV)**

Den här servicehandboken finns bara tillgänglig på engelska.

- Om en kunds servicetekniker har behov av ett annat språk än engelska, ansvarar kunden för att tillhandahålla översättningstjänster.
- Försök inte utföra service på utrustningen om du inte har läst och förstår den här servicehandboken.
- Om du inte tar hänsyn till den här varningen kan det resultera i skador på serviceteknikern, operatören eller patienten till följd av elektriska stötar, mekaniska faror eller andra faror.

**OPOZORILO
(SL)**

Ta servisni priročnik je na voljo samo v angleškem jeziku.

- Če ponudnik storitve stranke potrebuje priročnik v drugem jeziku, mora stranka zagotoviti prevod.
- Ne poskušajte servisirati opreme, če tega priročnika niste v celoti prebrali in razumeli.
- Če tega opozorila ne upoštevate, se lahko zaradi električnega udara, mehanskih ali drugih nevarnosti poškoduje ponudnik storitev, operater ali bolnik.

**DİKKAT
(TR)**

Bu servis kılavuzunun sadece ingilizcesi mevcuttur.

- Eğer müşteri teknisyeni bu kılavuzu ingilizce dışında bir başka lisandan talep ederse, bunu tercüme ettirmek müşteriye düşer.
- Servis kılavuzunu okuyup anlamadan ekipmanlara müdahale etmeyiniz.
- Bu uyarıya uyulmaması, elektrik, mekanik veya diğer tehlikelerden dolayı teknisyen, operatör veya hastanın yaralanmasına yol açabilir.

This page left intentionally blank.

Revision History

Rev	Date	Description
1	3-31-11	Initial release.

This page left intentionally blank.

Table of Contents

Chapter 1 INTRODUCTION.....	15
1 Introduction (Read First).....	15
1.1 Product Overview.....	15
1.2 Safety.....	16
1.3 Damage In Transportation.....	17
1.4 Product Locator.....	17
1.5 Installation and Setup Steering Guide.....	18
1.5.1 Scanner Interface Prerequisites.....	18
1.5.2 OR-Compatible Table Process.....	18
Chapter 2 INSTALLATION.....	21
1 Coil Cable Clamp Adapting Base Installation.....	21
1.4.1 Check M0032SS Coil Accessory Kit Contents.....	21
1.4.2 Locate Maquet Transfer Board.....	22
1.4.3 Installation.....	22
Chapter 3 ADJUSTMENTS AND CALIBRATIONS.....	25
1 Removing Cradle Assembly.....	25
2 Maquet Table Docking Pin Setup.....	27
Chapter 4 FUNCTIONAL CHECKS.....	31
1 Cradle Motion Test.....	31
2 Patient Table Movement Verification.....	34
Chapter 5 PARTS.....	37
1 Table and Cradle Assembly Spare Parts.....	37
Chapter 6 PLANNED MAINTENANCE.....	39
1 Planned Maintenance.....	39
Chapter 7 REPLACEMENTS.....	41
1 Safety Decal Location and Replacement.....	41

2 OR-Compatible Patient Transport Table Components.....	48
3 Cradle Wheel and Housing Replacement.....	49
4 Roller Blade Access for Wheel Cleaning or Replacement.....	52
5 Release Handle or Lock Lever Replacement.....	55
6 Surgical Table Dock Replacement.....	60
6.3.1 Dock Assembly Removal.....	61
6.3.2 Replacement Dock Assembly Installation.....	64
7 Touch and Go Strip and Side Bumper Replacement.....	66
Chapter 8 TROUBLESHOOTING.....	73
1 Table Down Inhibit Not Working When Docked.....	73

Chapter 1 Introduction

1 Introduction (Read First)

1.1 Product Overview

The OR-Compatible Table is a purchased option that provides the customer the ability to acquire image sets with the patient situated in a fixed surgical position. This OR table is used for Surgical Suite imaging only, and cannot be used for system calibrations because the flat surface cradle is not compatible with the calibration fixtures and phantoms. The system must be fully calibrated, using the primary Patient Transport Table, prior to the installation and setup of the OR-Compatible Table.

Illustration 1-1: MR750/MR450w OR-Compatible Table and Transfer Cradle



1.2 Safety



CAUTION

Potential for injury to service personnel!

The cradle used on the MR750/MR450w OR-Compatible Table weighs approximately 84 lb. (38 kg).

Use your judgment, as a second person may be required to assist in the cradle positioning process.



NOTICE

Follow ALL required safety and PPE procedures customary for your organization when working on this product.

- Wear gloves for hand protection.
- Wear safety glasses for eye protection.

1.3 Damage In Transportation

All packages should be closely examined at time of delivery. If damage is apparent, have notation "Damage in Shipment" written on all copies of the freight or express bill **before** delivery is accepted or signed for by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **must** be reported to the carrier immediately upon discovery or within **14 days** after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this 14-day period.

File a report by calling 1-800-548-3366, using Option 6. Contact your local Service Coordinator for more information on this process.

1.4 Product Locator

The Global Install Base Database (also known as Product Locator System) tracks shipment, trans-shipment, and field location of the serialized models. There are now two methods for submitting Product Locator information.

1. At this time, the preferred method of submitting information is the FE Site Verification website at: <http://gib.gehealthcare.com>

The FE Site Verification consists of three components available from the main menu:

- Install/deinstall product locator model and serial numbers
- Add/modify ship to address information
- Update CARES FE data for primary/secondary FEs

To obtain a copy of the FE training tutorial for using this website, a downloadable copy is available from the Product Locator Support Central webpage at: <http://supportcentral.ge.com/15563>

NOTE: One Shipping Card is filled out and submitted when shipped. (Extra cards are supplied for trans-shipments between storage and distribution points.) The installation card and extra shipment cards are attached.

2. Verify that serial and model numbers on each rating plate match the installation card numbers before removing installation cards. Process just the installation card and discard any extra shipment cards and labels.

NOTE: There may be one or more shipment cards and bar code labels with the installation card. These shipment cards are used to trace the transfer of serialized units between various inventory storage and distribution points until the product reaches its final installation destination.

3. Enter the information from the OR-Compatible Table installation card into the GIB database for this customer.

1.5 Installation and Setup Steering Guide

The OR-Compatible Table utilizes the same hardware interface to the MRI scanner for table docking and cradle attachment to the LPCA, so the majority of the service procedures will be the same. This document will only cover procedures specific to the OR-Compatible Table features, and provides assistance on installation, setup and service maintenance of the table.

1.5.1 Scanner Interface Prerequisites

Make sure the scanner interface hardware and adjustments are set up using the customer's primary Patient Transport Table. Refer to the Service Methods DVD that was shipped with the system. Perform all adjustments in the order listed below:

1. Leveling the Patient Transport
2. Table Top Height Adjustment
3. Cradle Release from Patient Transport Adjustment
4. Cradle Release Block Adjustment
5. Cradle Release from Carriage Adjustment
6. Cradle Release Pin Height Adjustment
7. Cradle Home Table Down Interlock Adjustment

1.5.2 OR-Compatible Table Process

For procedures marked with an asterisk (*), refer to the Service Methods DVD that was shipped with the system.

1. Perform these setup procedures on the OR-Compatible Table:
 - a. *Leveling the Patient Transport
 - b. *Table Top Height Adjustment
 - c. Removing Cradle Assembly
 - d. *Dual Table Alignment Procedure
 - e. *Cradle Release from Patient Transport Adjustment
 - f. *Cradle Release Block Adjustment
 - g. *Cradle Release from Carriage Adjustment
 - h. *Cradle Release Pin Height Adjustment

- i. *Cradle Home Table Down Interlock Adjustment
- j. Maquet Table Docking Pin Setup
- 2. Perform the following Functional Checks:
 - Cradle Motion Test
 - Patient Table Movement Verification

This page left intentionally blank.

Chapter 2 Installation

1 Coil Cable Clamp Adapting Base Installation

1.1 Personnel Requirements

Personnel Requirements	Preliminary Reqs	Procedure	Finalization
1	30 mins	30 mins	Not Applicable

1.2 Overview

This procedure is applicable to the Maquet Table Transfer Board which is used on all of the Maquet Surgical Tables or Maquet Transmobile Tables.

Commercial Catalog, M0032SS, Coil Accessory Kit, ships with every OR-Compatible Patient Transport Table Option order.

The purpose of this kit is to provide a means to secure the 6-Channel Flex Coil RF transmission leads to the bed of the Maquet Transfer Board. This is done to prevent injury to an immobilized patient if the coil's transmission leads were to get pulled during the diagnostic imaging process.

It is the responsibility of the GEHC Field Engineer to install the left and right hand adapting base assemblies onto the Maquet Transfer Board.

1.3 Preliminary Requirements

1.3.1 Tools and Test Equipment

Item	Qty	Effectivity	Part#	Manufacturer
Screwdriver; flat blade	1	-	5109895-10	-
Approved Utility knife	1	-	Local Purchase	-

1.3.2 Safety



NOTICE

Follow ALL required safety and PPE procedures customary for your organization when working on this product.

- Wear gloves for hand protection and safety shoes for foot protection
- Wear safety glasses for eye protection

1.4 Procedure

1.4.1 Check M0032SS Coil Accessory Kit Contents

1. Locate the package containing commercial catalog M0032SS.

2. Make sure there are no damaged or missing parts. Refer to the table below for content description and quantities.

Table 2-1: M0032SS Coil Accessory Kit Contents

Item	Qty	Part Number	Description
1	2	5260610	Coil Cable Clamp Assembly
2	1	5264719	Left Hand Adapting Base Assembly
3	1	5264719-2	Right Hand Adapting Base Assembly
4	2	5269784	Coil Cable Narrow Extension Assembly
5	10	5135972-19	M5x10, Screw, Metric, Nylon, Slotted Flat Head, ISO 2009, DIN 936 A

1.4.2 Locate Maquet Transfer Board


The Maquet Transfer Board is available to the customer in the following commercial catalog: M0071SS.

1. Locate the Maquet Transfer Board, Model 1180.71A0, which is shipped in catalog M0071SS.
2. Report any damaged or missing parts as required.

1.4.3 Installation

1. Position the Maquet Transfer Board, finished side down, on a clean surface. See below.


Illustration 2-1: Maquet Transfer Board Assembly

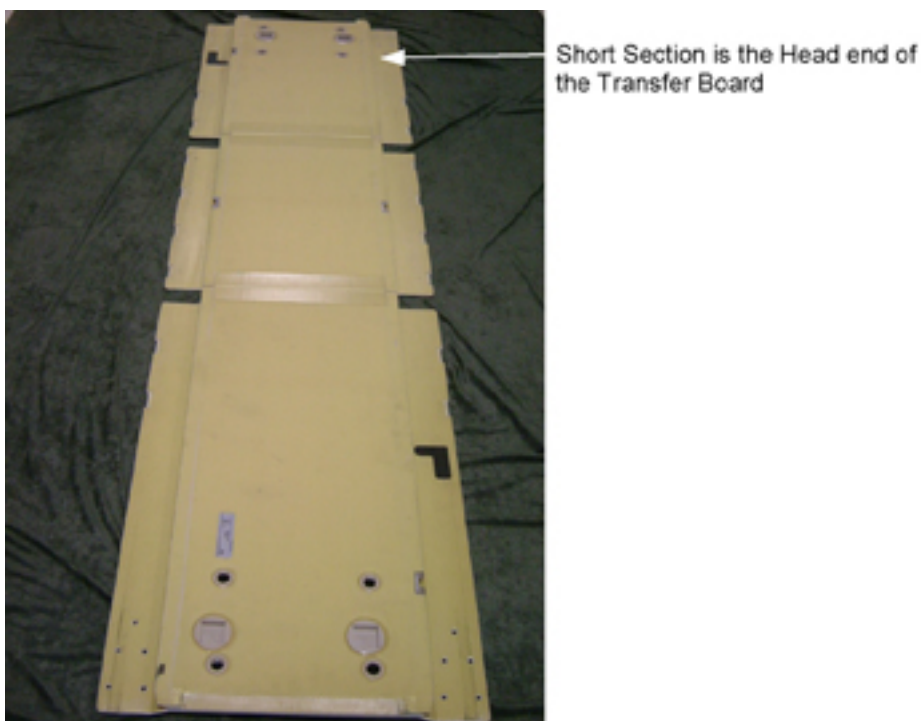

2. Mount the Left Hand Adapting Base Assembly, p/n 5264719, to the top surface of the Transfer Board using 5 plastic screws, p/n 5135972-19. See below.


NOTE: Three of the five screws will extend beyond the top surface of the Adapting Base Assembly due to the long length.

Illustration 2-2: Left Hand Adapting Base Assembly

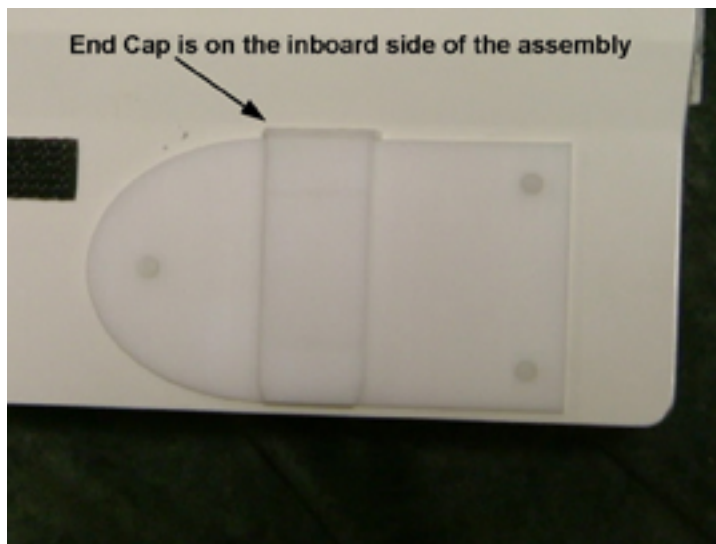

3. Make sure all screw head are tightly seated in the countersunk hole. See below.

Illustration 2-3: Transfer Board Bottom View

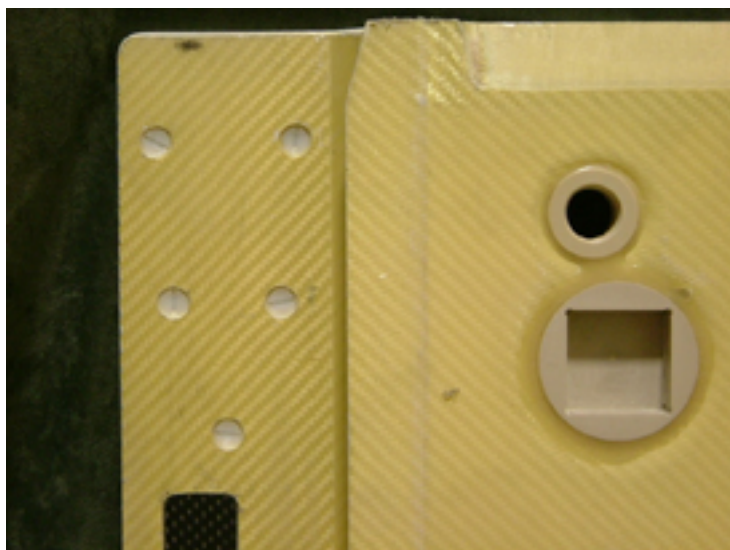

4. Trim off the excess screw shank flush with the adapting base using an approved utility knife. See below.

Illustration 2-4: Trim Excess Flush Using Approved Hand Protection and Cutting Tool



1.5 Finalization

Chapter 3 Adjustments and Calibrations

1 Removing Cradle Assembly

1.1 Personnel Requirements

Personnel Requirements	Preliminary Reqs	Procedure	Finalization
2		15 mins	


1.2 Overview

Removing the cradle assembly from the Patient Table requires at least two people to perform the task due to the weight (84 lbs) of the assembly.

1.3 Procedure

1. Undock the patient table.
2. Lock the patient table in place. Verify that the wheels are locked.
3. Person 1 and Person 2 should perform these tasks simultaneously.
 - a. Person 1: Manually rotate the transfer board cradle release lever (A).
 - b. Person 2: At the same time, press down on the cradle latch block (B) and push in the cradle release button (C) at the front of the cradle. Then raise the cradle about 3 to 4 inches. See illustration below.

Illustration 3-1: Cradle Release Mechanisms

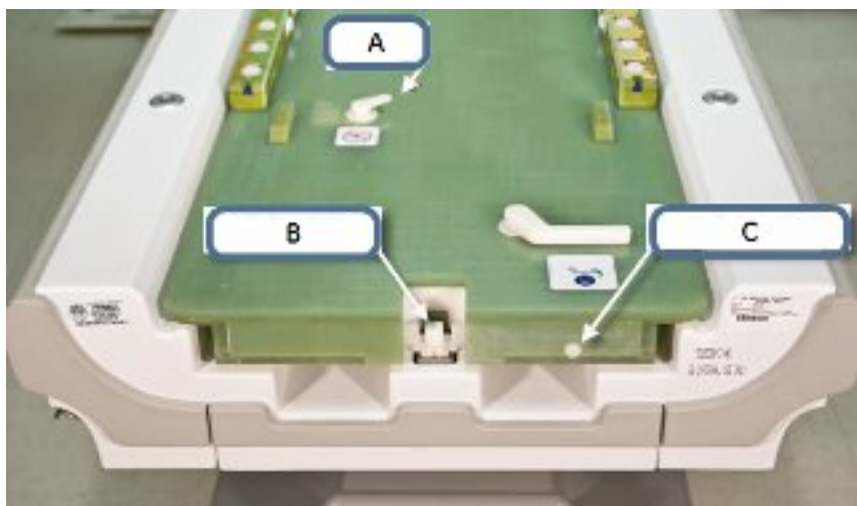

4. While doing all the tasks in step 3, slide the cradle forward enough so that the cradle covers the top area of the cradle latch block.

5. Person 2: Release the block and release the cradle button.

Illustration 3-2: Person 2



CAUTION

PINCHED FINGERS

Forward movement of the cradle can pinch your finger while holding down cradle latch block.

Move the cradle forward slowly and deliberately until the bottom edge of the cradle sufficiently covers the top of the cradle latch block, where you can remove your finger from the block.

6. Person 1 and Person 2 should perform these tasks simultaneously.
 - a. Person 1: Standing at the head end of the table, lift the cradle assembly off the table.
 - b. Person 2: Standing at the foot end of the table, lift the cradle assembly off the table.
7. Move the cradle to an appropriate area for servicing.

1.4 Finalization

No finalization steps.

2 Maquet Table Docking Pin Setup

2.1 Personnel Requirements

Personnel Requirements	Preliminary Reqs	Procedure	Finalization
1	Not Applicable	15 mins	Not Applicable

2.2 Overview

This procedure is applicable to the OR-Compatible Patient Transport Table Option. It requires one GE field engineer and one Maquet field engineer. The Maquet Surgical Table must have its height adjusted and set for patient transfer by a Maquet engineer prior to performing this procedure.

There is a docking pin located on the Maquet table that requires setup so it will properly activate a vertical locking mechanism within the OR-Compatible Table Cradle.


A brief description of the functionality of the table and cradle interlocks follows:

The Right and Left Vertical Lockers allow or restrict the movement of the Maquet Transfer Board.

- Right Vertical Locker - Activated by docking the OR-Compatible Table to a Maquet Surgical Table or Maquet Transmobile Table. A docking pin on the Maquet Table (s) presses a button on the cradle, which causes the Right Vertical Locker to retract. See below.
- Left Vertical Locker - Activated by Manual Release lever. When the user pulls the Manual Release lever it forces the vertical locker to retract. This action unlocks the Transfer Board from the cradle, which allows the Transfer Board to be moved onto the Maquet Surgical Table or Maquet Transmobile Table.

NOTE: When engaged, the Left and Right Vertical Lockers secure the Maquet Transfer Board to the OR-Compatible Cradle. Both Vertical Lockers must be retracted to allow the Transfer Board to be moved off of the cradle and onto a Maquet Table.

Illustration 3-3: Vertical Lockers In OR-Compatible Table Cradle

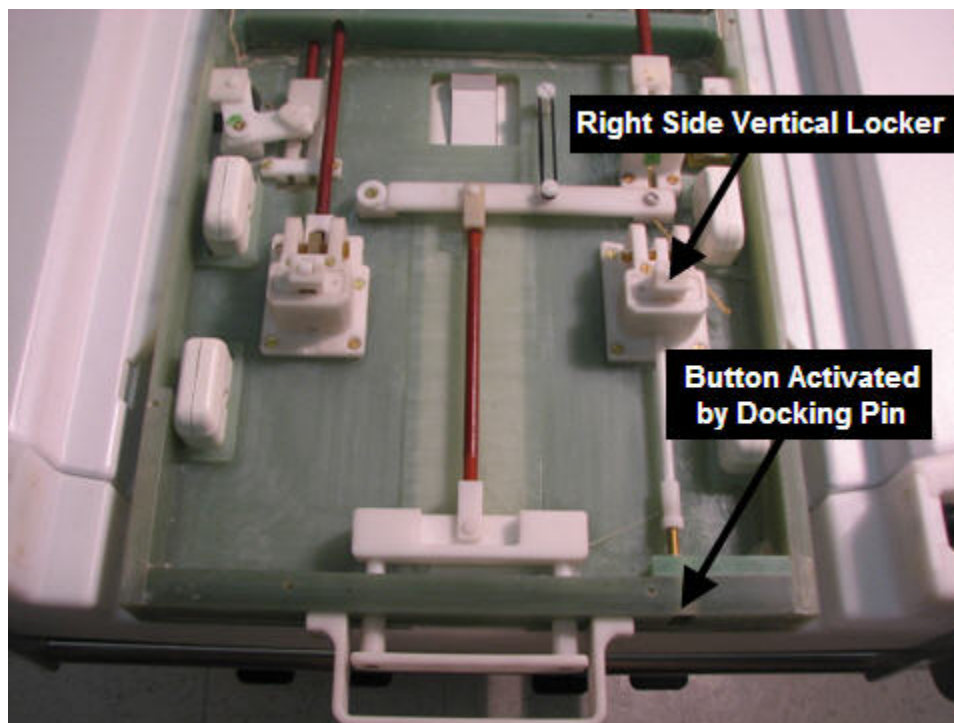


2.3 Preliminary Requirements

2.3.1 Tools and Test Equipment

Item	Qty	Effectivity	Part#	Manufacturer
Allen Key; 2.5 mm (part of tool set 5112581 or equivalent)	1	-	5109893-3	-

2.3.2 Safety

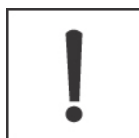

NOTICE

Follow ALL required safety and PPE procedures customary for your organization when working on this product.

- Wear gloves for hand protection
- Wear safety glasses for eye protection


2.4 Procedure

1. Move the OR-Compatible Table with Cradle assembly into the Operating Room and position it near to the Maquet Table but do not dock it.
2. With the assistance of the Maquet field engineer, set the Maquet table to the patient transfer height.

The next 2-steps are performed on the left front side of the Maquet Table.

3. Go to the left front side of the Maquet Table and locate the Docking Pin Assembly. See below.

Illustration 3-4: Maquet Table Docking Pin Assembly

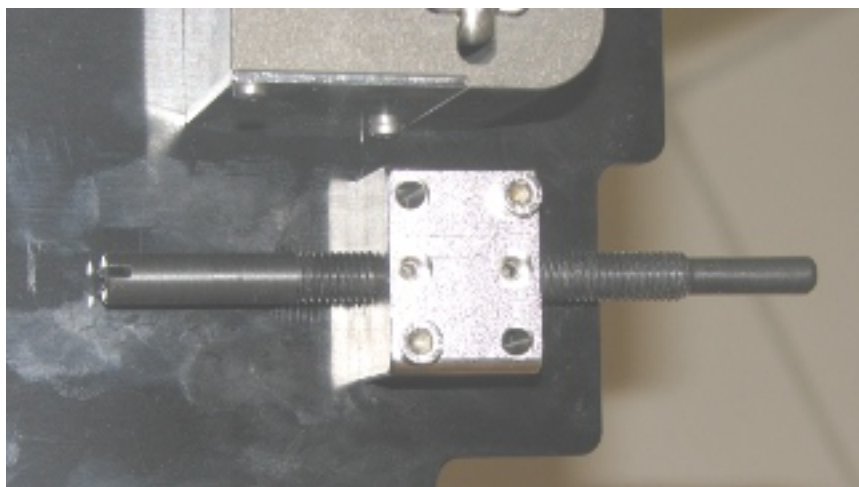

4. Loosen the Docking Pin Set Screws, but do not remove them. Verify the docking pin can be easily moved within the mounting assembly.
5. Move the OR-Compatible Table, with Cradle, in front of the Maquet Table and dock it.

NOTE: The table dock mechanism will automatically engage and lock the table into place.

Illustration 3-5: Table Lock Mechanism



6. Manually push the docking pin forward so it engages with the button on the GE Cradle and continue to apply enough force to cause the Right Vertical Locker on the GE Cradle to retract.
7. Once retracted, tighten the setscrews to lock the Docking Pin in place.
8. Verify the Right Vertical Locker remains retracted.
9. Rotate the Lock Knob to unlock the OR-Compatible Table from the Maquet Table, Separate the two tables and verify the Right Vertical Locker has returned to the extended position.
10. Re-dock the two tables and verify the Right Vertical Locker moves to the retracted position.
11. Perform this procedure on any Maquet Transmobile Tables that will be used by the customer for this surgical suite.

2.5 Finalization

No finalization steps.

Chapter 4 Functional Checks

1 Cradle Motion Test

1.1 Personnel Requirements

Personnel Requirements	Preliminary Reqs	Procedure	Finalization
1	Not Applicable	30 mins	Not Applicable

1.2 Overview

This test moves the cradle into the magnet bore, and then returns it to its starting position. When this test is run, the cradle is automatically moved into the magnet bore a maximum distance of 600 mm, then is returned to its starting position. Any error messages generated during this test are stored in the Error Message Log. If the cradle starts in the Home position, it is moved about 100 mm into the magnet bore, is paused there, then is advanced a maximum of 600 mm into the bore, then is returned close to the 100 mm position. If the cradle starts in a position beyond 1700 mm into the bore (max. Is 1160 for LCC), it is automatically returned to the Home position, then is advanced a maximum of 600 mm into the bore, then is returned close to the 100 mm position.

1.3 Preliminary Requirements

1.3.1 Safety



⚠ DANGER

STRONG MAGNETIC FIELD!

FERROUS MATERIALS CAN BECOME DANGEROUS PROJECTILES IN THE PRESENCE OF THE MAGNETIC FIELD PRODUCED BY THE SIGNA MAGNET. DO NOT BRING ANY FERROMAGNETIC TOOLS OR EQUIPMENT INTO THE MAGNET ROOM.



⚠ CAUTION

POSSIBLE PERSONAL INJURY AND/OR EQUIPMENT DAMAGE.


OBSTRUCTIONS OR PERSONNEL IN CRADLE PATH MAY CAUSE PERSONAL INJURY OR EQUIPMENT DAMAGE.

VERIFY THAT MAGNET BORE IS CLEAR OF PERSONNEL AND EQUIPMENT BEFORE STARTING CRADLE MOTION TEST.

1.4 Procedure

1. Start the Service Browser and go to the MR Service Desktop on the host computer.
2. Go to the [Diagnostics tab], [Peripherals], [SRI Functional Tests], and select the [Motion Test]. See below.

Illustration 4-1: SRI Functional Tests

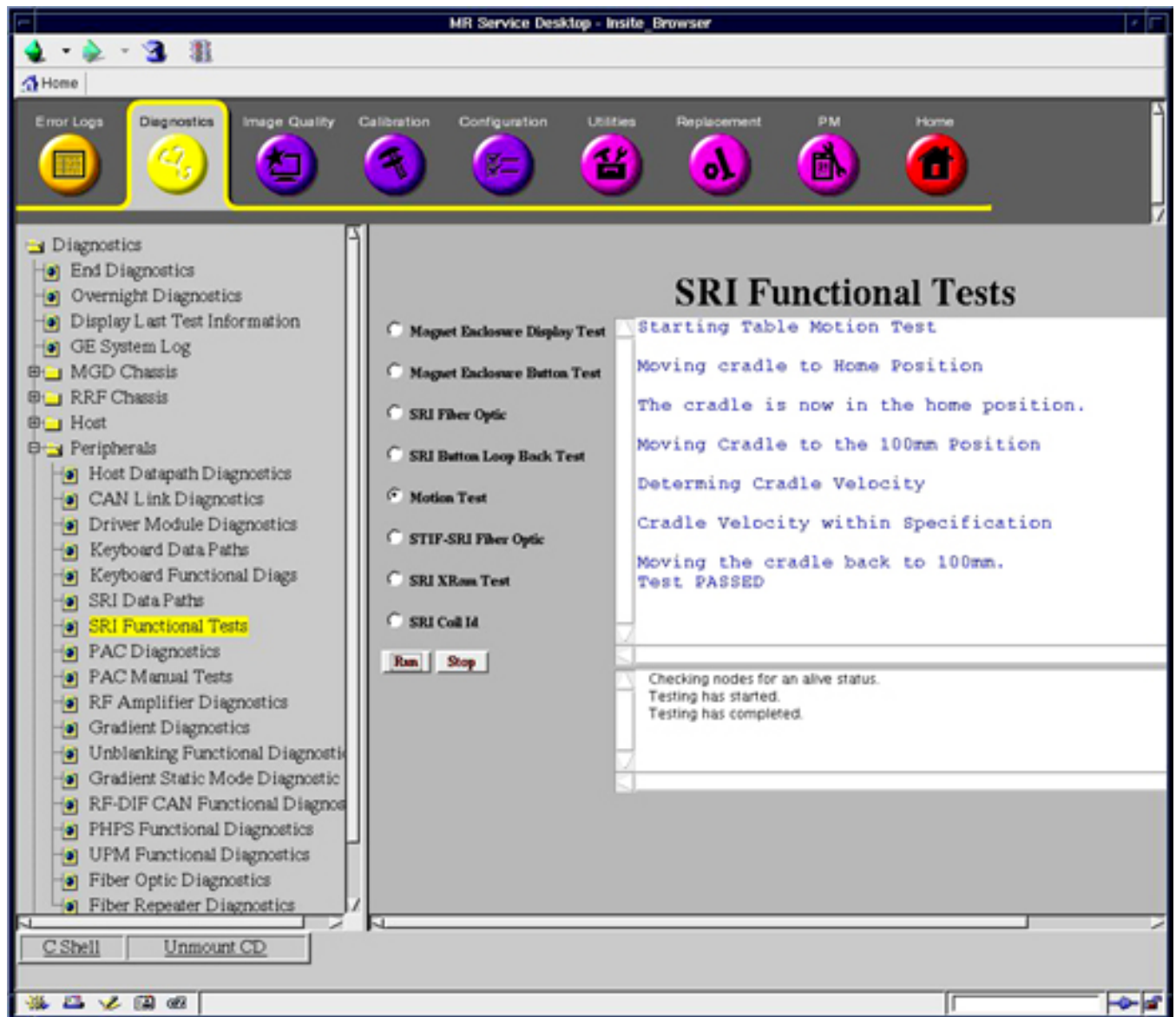

3. Select [Motion Test] and then click [Run].
4. A pop-up window will appear warning that the cradle is about to move. Press [Yes] to continue.
5. Verify that the cradle moves smoothly, and that magnet enclosure cradle position display increments, or decrements, corresponding to cradle travel into or out of magnet bore.
6. After the cradle stops, verify the test passed in the status window of the diagnostic tool.
7. Go to [Diagnostics], End Diagnostics and select the [End Diagnostics] Button to exit the diagnostics mode.

1.5 Finalization

No finalization steps.

2 Patient Table Movement Verification

2.1 Personnel Requirements

Personnel Requirements	Preliminary Reqs	Procedure	Finalization
1	Not Applicable	20 mins	Not Applicable

2.2 Overview

The OR-Compatible Patient Table may be operated in either the docked or undocked position. Undocked, it may be raised or lowered via the foot pedals located at the rear of the table. In addition, floor locks are available for locking the wheels during patient transfer. While the OR-Compatible Patient Table is docked to the magnet enclosure, the table may not be driven up or down by either the dock pedals or the foot pedals at the rear of the table.

2.3 Procedure

1. Undock patient table by pressing second pedal from left at rear of table.
2. Press the Table Down pedal (far right) to initiate downward movement of the table top. Downward movement of table top should stop when pedal is released. Verify that table travels from highest to lowest position in no more than 15 seconds with table unloaded.
3. Pump the Table Up pedal (second from right) to initiate upward movement of table. Verify that the table reaches its highest position in 15 pedal strokes or fewer. For accurate check of this function to specification given, pedal strokes must be complete (i.e. push pedal all the way down and let pedal come all the way up before starting next downward movement).
4. Press the Caster Lock pedal on one of the rear casters to lock both rear casters. See below.

Illustration 4-2: Engaging the Caster Lock Pedal




- a. Verify caster locks inhibit all Patient Table movement when pedal is in locked position.
- b. Verify there is no movement of cradle while Patient Table is undocked.

**NOTICE**

The OR-Compatible Table must always be in its highest position when docking to the scanner.


5. Dock Patient Table by aligning table with dock and then pressing the far left dock pedal (see below). Successful docking should not require excessive force. When docked, verify that carriage advances forward until it connects with cradle latch immediately after docking.

Illustration 4-3: Dock Pedal on Patient Table

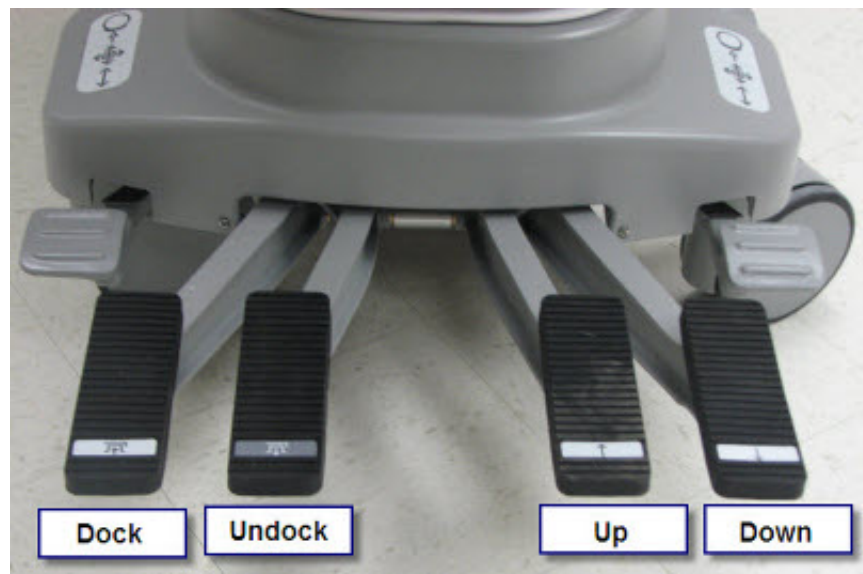

6. Verify that the emergency release lever (located on the right side of the dock) releases the Patient Table from the dock when it is pulled (make sure LPCA has been disengaged from cradle). See below.

Illustration 4-4: Dock Emergency Release Handle

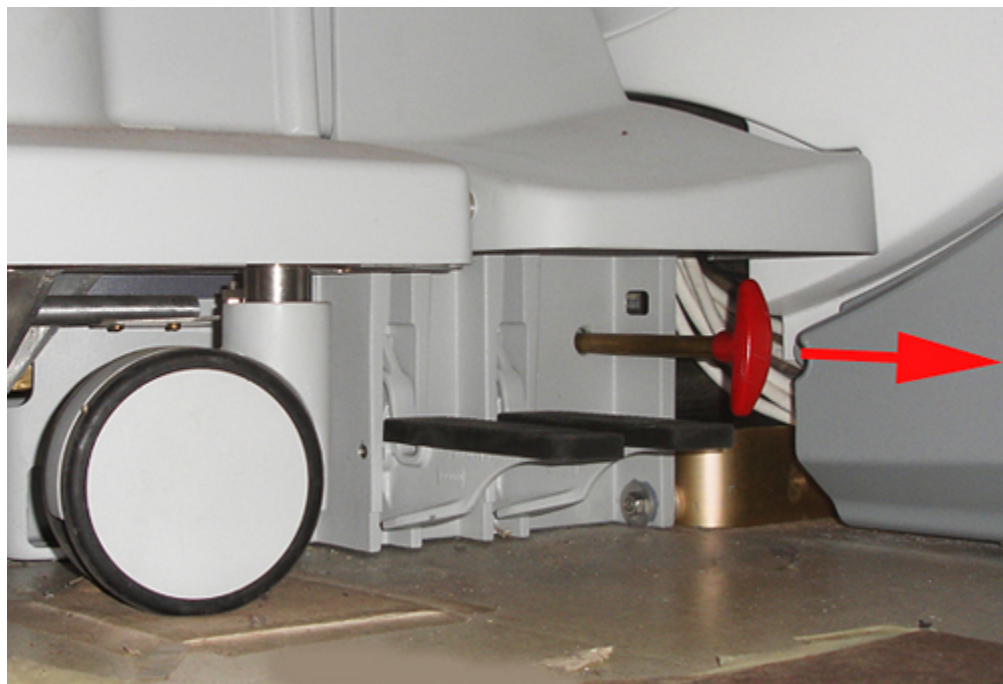

7. Redock the patient table.
8. Verify that when you squeeze the cradle release handle at the rear of the cradle, it disengages cradle from LPCA. Refer to the Service Methods DVD shipped with the system for procedure for Cradle Release Adjustments (Cradle Release Block Adjustment, Cradle Release From Carriage Adjustment, or Cradle Release from Patient Table).
9. With the cradle fully in the bore, verify that the Table Down pedals on both the dock and table itself do not allow the table to go down. If the table does move down, perform Cradle Home Table Down Interlock Adjustment (Refer to the Service Methods DVD shipped with the system for this procedure).
10. Move the cradle to the home position. Verify that the Table Down pedals on both the dock and table itself do not allow the table to go down. If the table does not move down, refer to Troubleshooting (Table Down Inhibit Mechanism) in this manual.

2.4 Finalization

No finalization steps.

Chapter 5 Parts


1 Table and Cradle Assembly Spare Parts

This procedure is applicable to the OR-Compatible Patient Transport Table Option.

Table 5-1: OR-Compatible Patient Transport Table New Spare Parts

Item	Name	Part Number	Where Used
1	Reference Template for checking compatibility of articles on the patient table to 60 cm clear bore opening.	5389349	MR750 OR-Compatible Table
2	Reference Template for checking compatibility of articles on the patient table to 70 cm clear bore opening.	5389349-2	MR450w OR-Compatible Table
3	OR-Compatible Cradle Assembly with packaging crate	5367845-2	OR-Compatible Table MR750/MR450w
4	Upper Table Dock Assembly with packaging crate	5363547-2	MR750/MR450w OR-Compatible Table
5	Mayfield Pad	5391776	For the surgical head clamp
6	OR-Compatible Table Decals	5390451	Used in various locations on the OR-Compatible Table and Dock
7	Cradle Decals	5371032	OR-Compatible Cradle Assembly
8	Maquet Roller Blade Assembly	5367850	Cradle Assembly
9	Stop and Lever Kit	5390460	Cradle Assembly
10	Maquet Roller	5367873	Cradle Assembly
11	Cradle Wheel Kit includes set of six wheel housings and wheels. Existing part used on Primary Patient Transport Table Cradle	5179918	Cradle Assembly
12	Rear Side Bumper LH	5371065	OR-Compatible Table MR750/MR450w
13	Rear Side Bumper RH	5371065-2	OR-Compatible Table MR750/MR450w

Illustration 5-1: OR-Compatible Patient Transport Table And Cradle Components

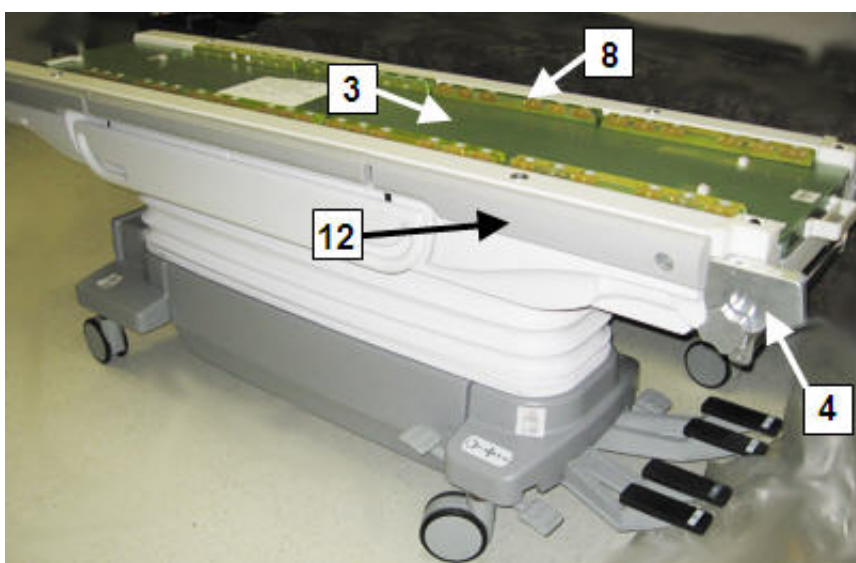

This page left intentionally blank.

Chapter 6 Planned Maintenance

1 Planned Maintenance

The OR-Compatible Table follows the same planned maintenance schedule as the Patient Transport Table used on the system.

Refer to the system DVD for schedule and tasks.

This page left intentionally blank.

Chapter 7 Replacements

1 Safety Decal Location and Replacement

1.1 Personnel Requirements

Personnel Requirements	Preliminary Reqs	Procedure	Finalization
1	Not Applicable	Not Applicable	Not Applicable

1.2 Overview

This procedure is applicable to the MR750/MR450w OR-Compatible Patient Transport Table Option.

The safety decals used on this table identify operation and heavy cradle weight. These should be replaced when they are found missing or unreadable.

Table Cradle Decals Locations

Illustration 7-1: Manual Operation Inhibit



Illustration 7-2: Transfer Board Release

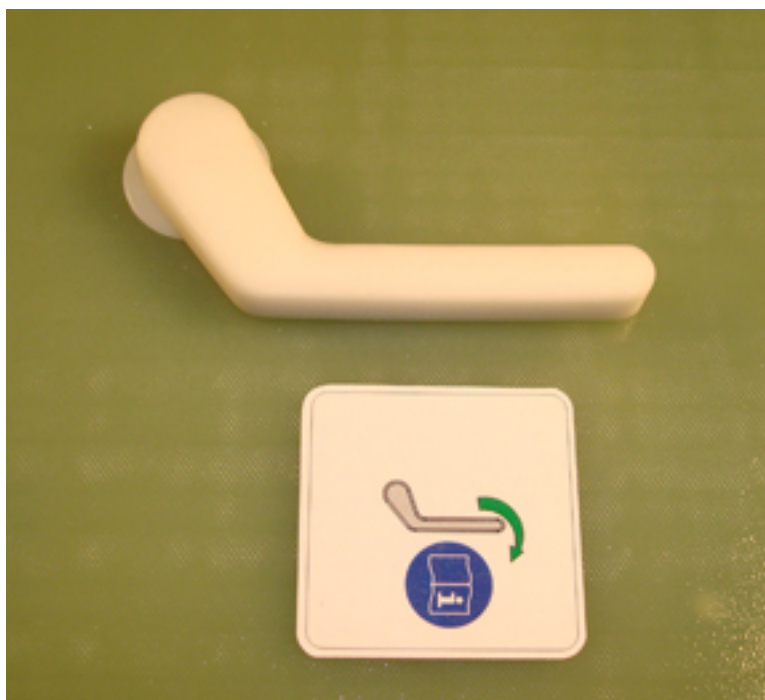


Illustration 7-3: Heavy Objects / Cradle Release

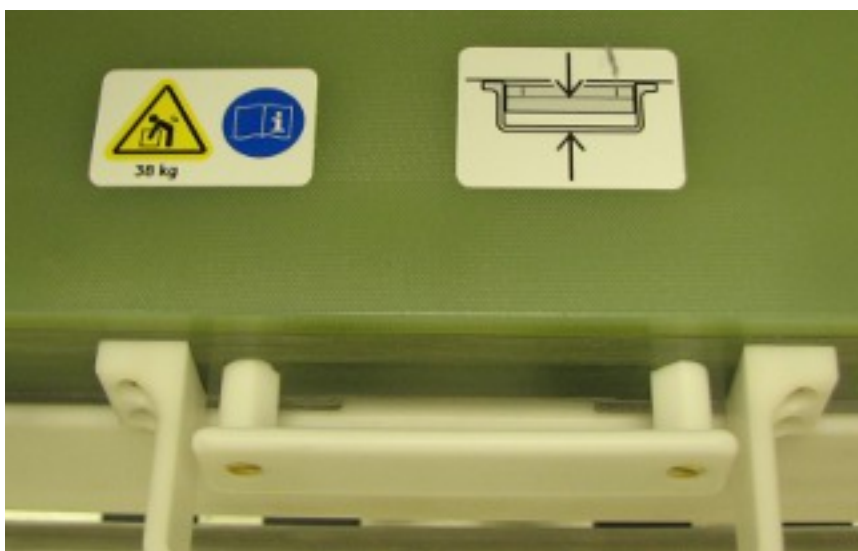


Table Decals Locations

Illustration 7-4: Maquet Dock Locking Mechanism



Illustration 7-5: Table Weight Warning

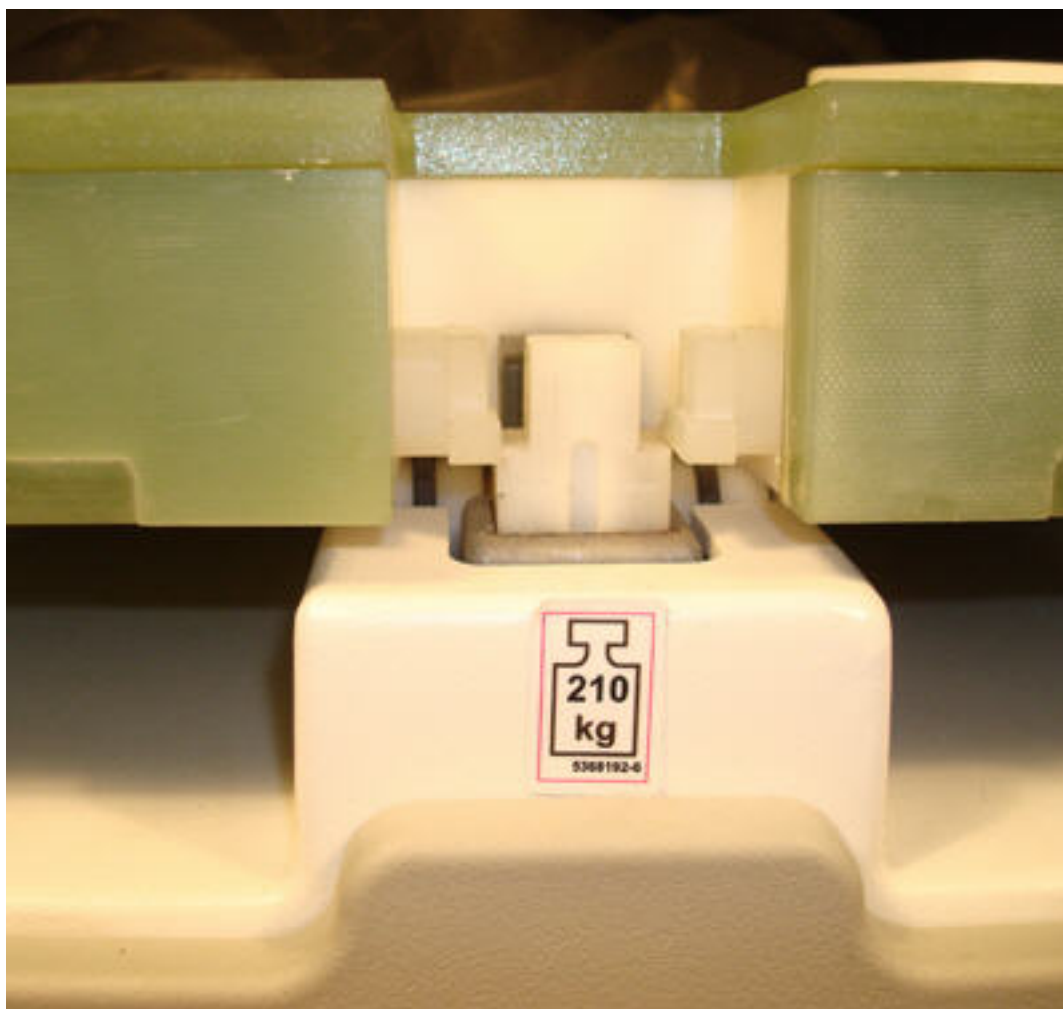


Illustration 7-6: Table Down Operation - Table Pedal End



Illustration 7-7: Table Down Operation - Scanner End



1.3 Preliminary Requirements

1.3.1 Consumables

Item	Qty	Effectivity	Part#	Manufacturer
Isopropyl Alcohol	As required	-	-	-
Towels or Cloths	As required	-	-	-

1.3.2 Replacement Parts

Item	Qty	Effectivity	Part#	Manufacturer
Cradle Decals	1	-	5371032	-
OR-Compatible Table Decals	1	-	5390451	-

1.3.3 Safety



NOTICE

Follow ALL required safety and PPE procedures customary for your organization when working on this product.

- Wear gloves for hand protection
- Wear safety glasses for eye protection

1.4 Procedure

1. Undock the Patient Transport Table from the magnet and move it to a location suitable for servicing.
2. Clean the area of interest using a clean towel and Isopropyl alcohol. Remove any pieces of the old decal. Allow the surface to dry before applying the new decal.
3. Remove the protective backing from the decal and install the decal. Make sure to smooth out any air bubbles.

1.5 Finalization


No finalization steps.

2 OR-Compatible Patient Transport Table Components

Several procedures relating to components of the Patient Table require the removal of the table Upper Bellows and Lower Base Side Cover. Refer to the Service Methods DVD for more information.

The following illustration depicts the major replaceable components for the MR750 / MR450w OR-Compatible Patient Table.

Illustration 7-8: Patient Table Main Components

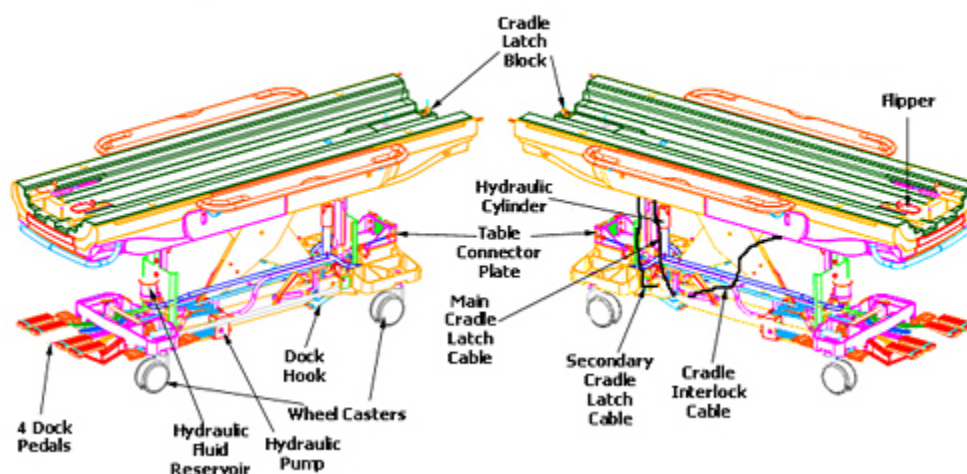

3 Cradle Wheel and Housing Replacement

3.1 Personnel Requirements

Personnel Requirements	Preliminary Reqs	Procedure	Finalization
1 or 2, depending on system	Not Applicable	30 mins	Not Applicable

3.2 Overview

This procedure is applicable for the MR750 / MR450w OR-Compatible Patient Transport Table.

The cradle used on the MR750 / MR450w OR-Compatible Patient Transport Table approximately 84 lbs. (38 kgs.), which may require 2 persons for placing the cradle into a serviceable position.

3.3 Preliminary Requirements

3.3.1 Tools and Test Equipment

Item	Qty	Effectivity	Part#	Manufacturer
Non-Magnetic Hand Tool Kit	1	-	5113258 or 5112581	-
Flat Blade Screwdriver (from kit)	1	-	-	-
Pliers (from kit)	1	-	-	-

3.3.2 Replacement Parts

Item	Qty	Effectivity	Part#	Manufacturer
Wheel and Housing Kit (containing 6 Wheel and Housing assemblies)	As required	-	5179918	-

3.3.3 Safety



CAUTION

POTENTIAL FOR INJURY TO SERVICE PERSONNEL.

THE CRADLES USED ON THE MR750 / MR450w OR-COMPATIBLE TABLE WEIGH APPROXIMATELY 84 LBS. (38 KGS.).

THE CRADLE CAN BE TURNED OVER AND SUPPORTED BY THE TABLE FOR THIS REPLACEMENT PROCEDURE. USE YOUR JUDGMENT AS A SECOND PERSON MAY BE REQUIRED TO ASSIST IN THE CRADLE POSITIONING PROCESS.

**NOTICE**


Follow ALL required safety and PPE procedures customary for your organization when working on this product.

- Wear gloves for hand protection.
- Wear safety glasses for eye protection.

3.4 Procedure

1. Undock the OR-Compatible Patient Transport Table from the magnet and move it to a suitable location for servicing.
2. Lock the table castors to prevent table movement during servicing.
3. Release the cradle from the Patient Transport Table. See [Chapter 3, Removing Cradle Assembly](#).
4. Flip over the cradle to access the Wheel and Housing assemblies underneath. See below.


Illustration 7-9: Remove the Wheel and Housing

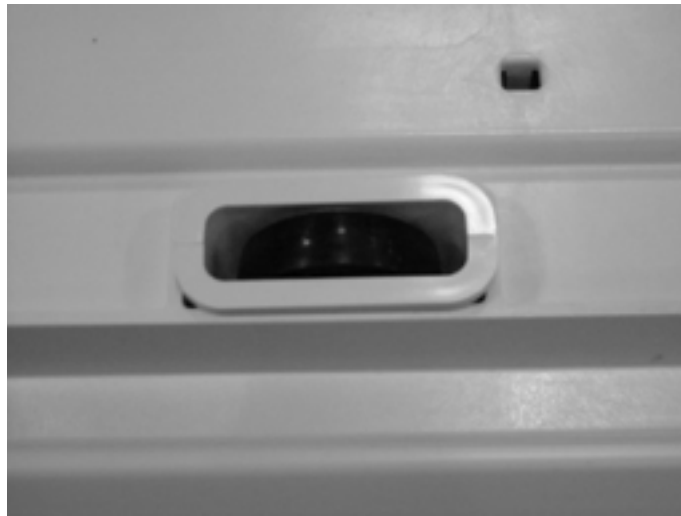

NOTE: The Wheel and Housing assembly is press fit into the pocket on the cradle. The choice of tools for removing the Wheel and Housing is dependent on how tight the fit is.

5. Pry and lift the defective Wheel and Housing assembly out from the pocket.
6. Install the new Wheel and Housing assembly.

Illustration 7-10: Wheel and Housing Assembly

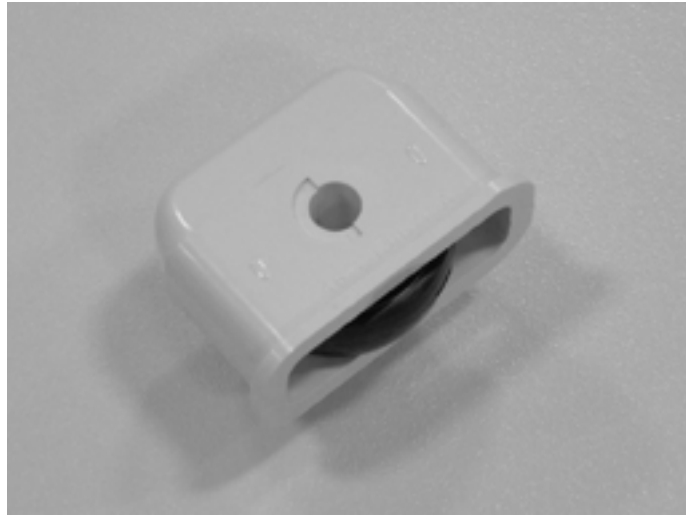

7. Verify the wheel spins freely within the housing.
8. Replace other defective Wheel and Housing assemblies as required.

3.5 Finalization

No finalization steps.

4 Roller Blade Access for Wheel Cleaning or Replacement

4.1 Personnel Requirements

Personnel Requirements	Preliminary Reqs	Procedure	Finalization
1	Not Applicable	1 hour	10 mins

4.2 Overview

This section is applicable to the OR-Compatible Patient Table cradle assembly. Follow it to access the wheels housed within the Roller Blade assembly for cleaning or replacement purposes.

4.3 Preliminary Requirements

4.3.1 Tools and Test Equipment

Item	Qty	Effectivity	Part#	Manufacturer
No. 2 Phillips Screwdriver (included in Titanium Tool Set)	1	-	5113258 or 5112581	-
Flat-Blade Screwdriver (included in Titanium Tool Set)	1	-	5113258 or 5112581	-
1/16 in. Allen Hex Key	1	-	Field Supplied	-

4.3.2 Consumables

Item	Qty	Effectivity	Part#	Manufacturer
Isopropyl Alcohol	1	-	Field Supplied	-
Cotton Swabs	20	-	Field Supplied	-
Paper Towels	1 roll	-	Field Supplied	-

4.3.3 Replacement Parts

Item	Qty	Effectivity	Part#	Manufacturer
Roller Blade Assembly	1	-	5367850	-
Maquet Roller	1	-	5367873	-

4.4 Procedure

1. Undock the table from the scanner and move to an appropriate area for servicing.
2. Lock the table castors to prevent table movement during servicing.

3. Determine which bank of roller blades requires servicing. See [Illustration 7-11](#).

NOTE: There are four banks of roller blade assemblies on each side of the cradle. Each assembly contains both large and small wheels.


4. Release the cradle and move it forward about 12 inches.
5. Tilt the cradle on its side to allow access to the appropriate roller blade assembly.

Illustration 7-11: Tilting Cradle to Access Blade Assembly



6. Hold down the roller blade release button with a flat-head screwdriver, and push the roller blade assembly away from the slotted pegs holding it to the cradle (shown above).
7. Remove the two screws holding the blade cover to the assembly to access the desired wheel.

Illustration 7-12: Accessing Marquet Roller

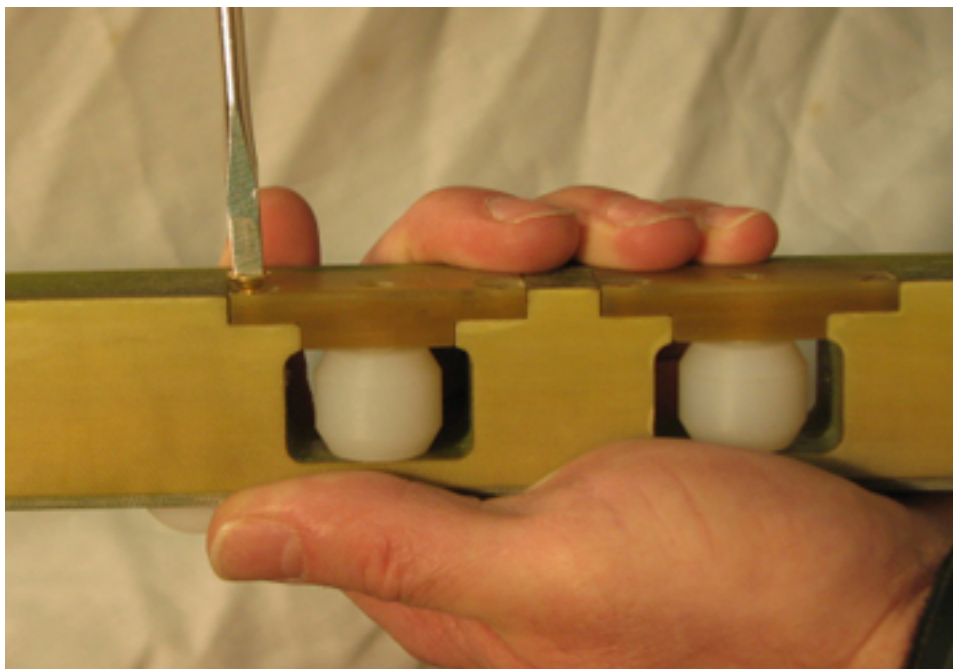

8. Use a cotton swab and Isopropyl alcohol to clean the wheel cavity and the hole that supports the wheel shaft.

NOTE: Do not use any lubricant in the shaft hole or on the wheels. All surfaces must be clean and dry before inserting wheels.

9. Reinstall the wheels and the assembly by following the reverse order of these instructions.

4.5 Finalization

1. Verify that all wheels rotate smoothly by running your hand over the top of them.
2. Return the cradle to its normal operating position once all cleaning or repair is completed.

5 Release Handle or Lock Lever Replacement


5.1 Personnel Requirements

Personnel Requirements	Preliminary Reqs	Procedure	Finalization
1	Not Applicable	Not Applicable	Not Applicable

5.2 Overview

The Release Handle and the Lock Lever are located at the head end of the Cradle. Each assembly has a seal and an o-ring to prevent fluids from leaking into the Cradle. See below.

Illustration 7-13: Release Handle and Lock Lever, End View

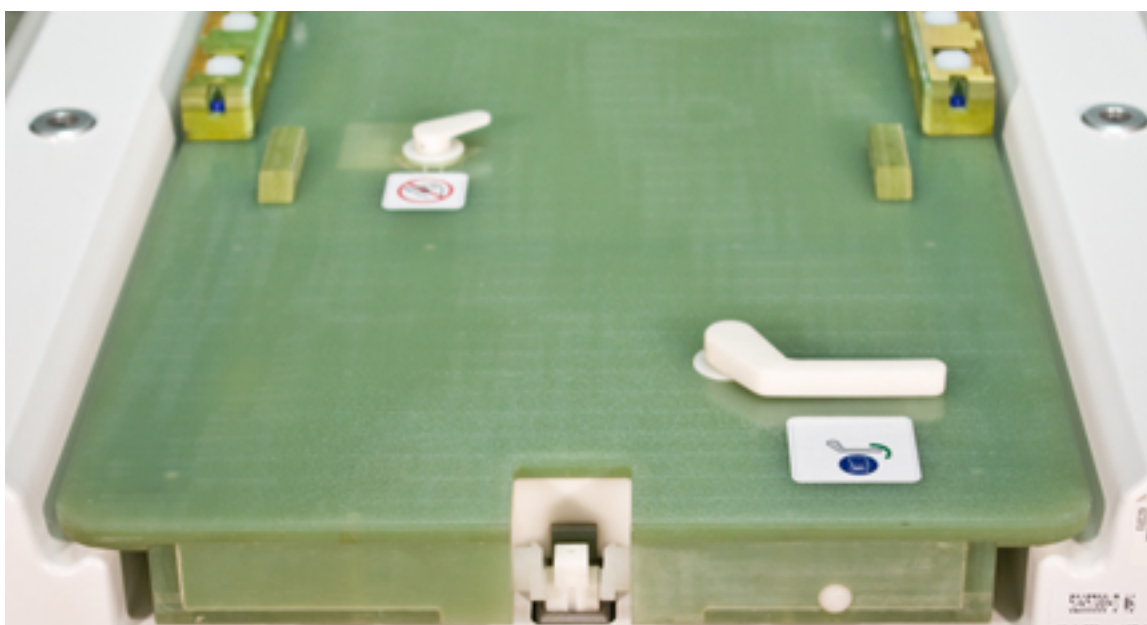

The Release Handle, operated by the customer, retracts a vertical locker that is engaged into the bottom of the Transfer Board. When this locker is retracted it allows the Transfer Board to be manually moved from the OR-Compatible Table to the Maquet Table.

The Lock Lever, operated by the Transfer Board, retracts a side locker that is engaged into a block on the topside of the table. When this locker is retracted it allows the Transfer Board and Cradle to be moved in and out of the magnet bore.

5.3 Preliminary Requirements

5.3.1 Tools and Test Equipment

Item	Qty	Effectivity	Part#	Manufacturer
Non-Magnetic Hand Tool Kit	1	-	5113258	-

Item	Qty	Effectivity	Part#	Manufacturer
Screwdriver, 1/8-inch, flat blade	1	-	-	-

5.3.2 Replacement Parts

Item	Qty	Effectivity	Part#	Manufacturer
Stop and Lever Kit (Kit contents listed below)	1	-	5390460	-
Manual Release Handle	1	-	5371205	-
Seal	2	-	5270414	-
Right Side Locker Lever	1	-	5376139	-
O-ring; Vitron Parker # 2-014; ID 12.42 mm; Parker Compound V0747-75	2	-	5304495	-
Screw Set 4-40 x .187 long; pure brass	2	-	5305074	-
Tube - RTV-162, Non-corrosive Adhesive	1	-	46-170137P1	-

5.3.3 Safety



NOTICE


Follow ALL required safety and PPE procedures customary for your organization when working on this product.

- Wear gloves for hand protection
- Wear safety glasses for eye protection

5.4 Procedure

1. Loosen the set screw located on the side of the Release Handle or Lock Lever. See below.

Illustration 7-14: Release Lever Assembly

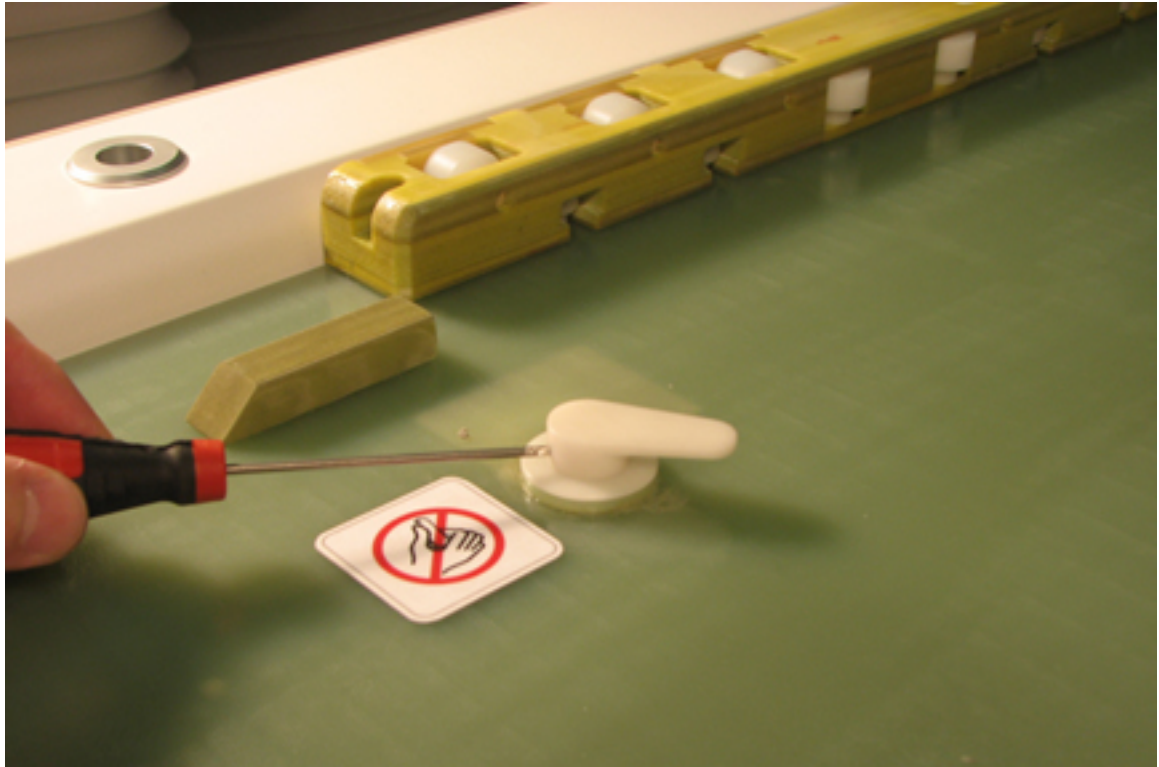


2. Remove the Release Handle or Lock Lever from the shaft.
3. Using a flat blade screwdriver, pry up the seal. See below.

Illustration 7-15: Pry Up the Seal



4. Remove the O-ring from the seal. See below.

Illustration 7-16: Remove the O-Ring from the Seal

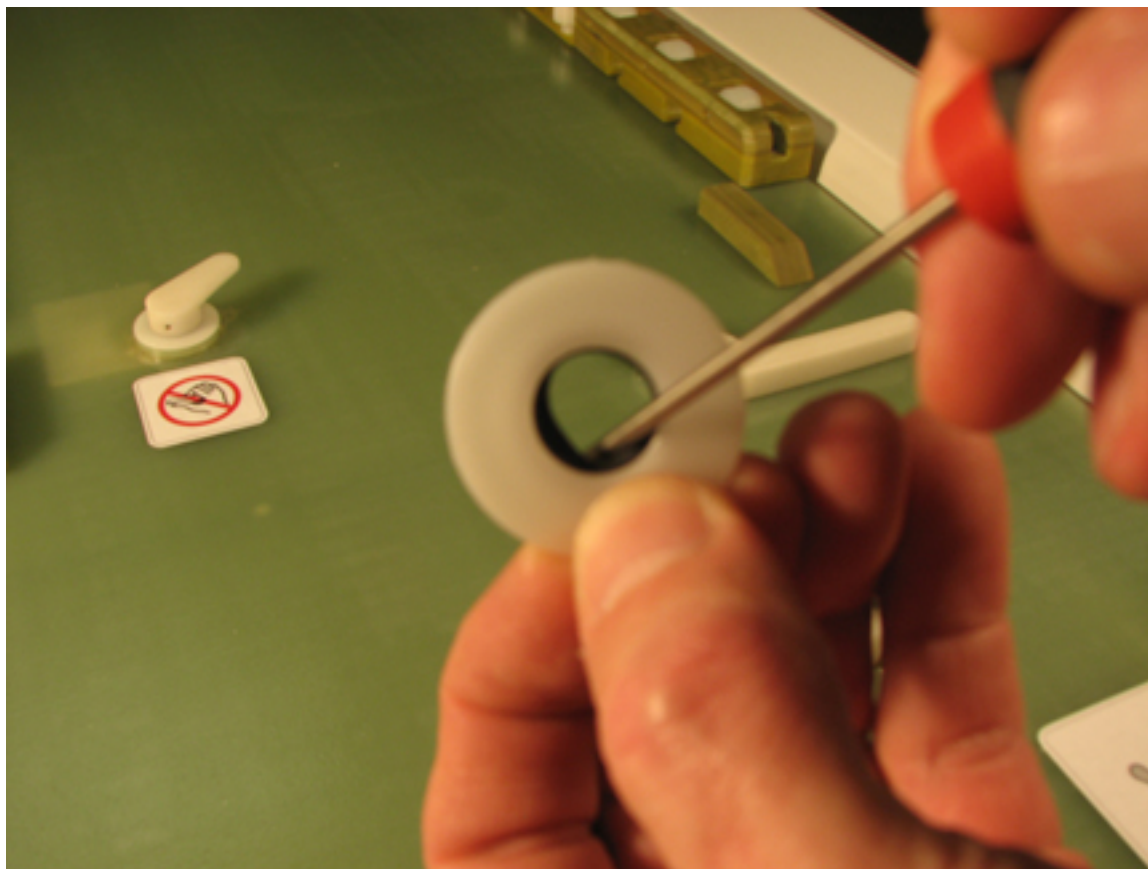

5. Install a new O-ring. See above.
6. Add a small amount of RTV around the seal and install around shaft.
7. Clean up any RTV that may have squeezed out during the seal installation.
8. Reinstall the Release Handle or Lock Lever and tighten the set screw.
9. Verify that the lever is secured to the shaft. Operate the lever to verify the operation of the locker at the other end of the cradle.

5.5 Finalization

No finalization steps.

6 Surgical Table Dock Replacement

6.1 Personnel Requirements

Personnel Requirements	Preliminary Reqs	Procedure	Finalization
1		90 mins	

6.2 Preliminary Requirements

6.2.1 Tools and Test Equipment

Item	Qty	Effectivity	Part#	Manufacturer
5/32 in. Allen Key	1	-	-	-
No. 2 Phillips-head Screwdriver	1	-	-	-
3/8 in. Combination Open Box Wrench	1	-	-	-

6.2.2 Consumables

Item	Qty	Effectivity	Part#	Manufacturer
2 in. wide Duct Tape	As needed	-	-	-

6.2.3 Replacement Parts

Item	Qty	Effectivity	Part#	Manufacturer
OR-Compatible Patient Table Dock Assembly	1	-	5363547	-

6.2.4 Safety



CAUTION

Exposure to Pathogens!

Potential pathogen fluids may be present on the table. If fluids are present, notify hospital personnel for cleanup and disposal.

Wear nitrile gloves and safety glasses while servicing the OR-Compatible Table.

6.2.5 Required Conditions


Condition	Reference	Effectivity
Access to the Operating Room is required to align the replaced table dock to the Operating Room Maquet Surgical Table.	-	-

6.3 Procedure

6.3.1 Dock Assembly Removal


1. Undock the OR-Compatible Table and move it out of the Magnet Room.
2. Set the brakes on the casters to prevent the table from moving.
3. Move the cradle forward approximately one foot to allow access to the dock. Refer to Removing Cradle Assembly for details.
4. Remove the Rear Bottom Cover by removing the four Phillips-head screws shown below.

Illustration 7-17: Rear Bottom Cover

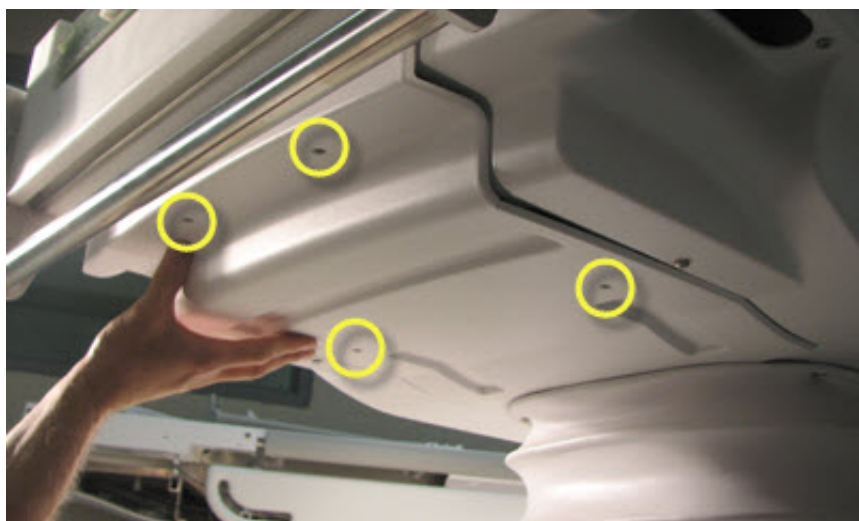


5. Remove the left and right Upper Rear Side Covers. (There are seven Phillips-head screws on each cover.)
6. Remove the left and right Rear Side Bumpers, secured with 3/8-inch hex nuts and brass washers, and the Rear Bottom Cover.

Illustration 7-18: Table Covers and Bumper

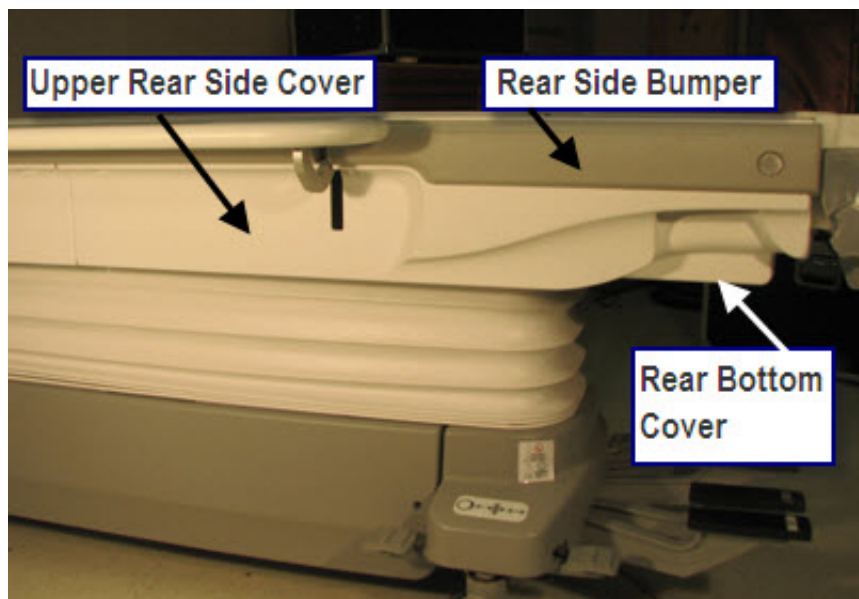

7. Secure the Upper Table Dock to the bezel by wrapping duct tape around the end of the table and under the Table Dock Bar.

Illustration 7-19: Table Dock Secured with Tape

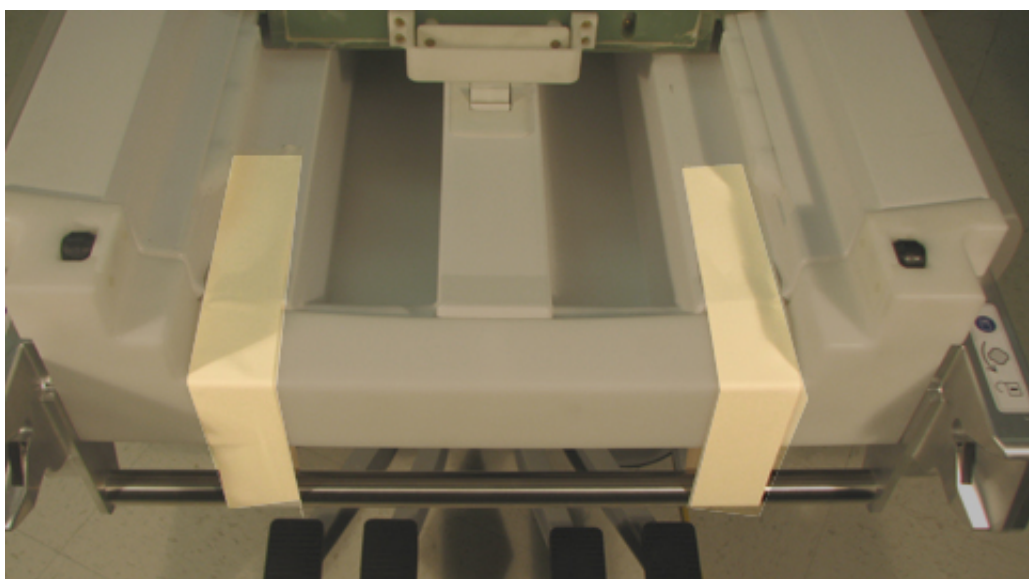


8. Remove the three Allen-head screws from the Upper Table Dock via the access holes on each side of the table.

Illustration 7-20: Removing Allen Screws from Dock



9. While supporting the dock assembly with one hand, remove the strips of tape securing the dock to the table bezel.
10. Lower the dock assembly and remove it from the table as shown below.


Illustration 7-21: Lowering and Removing Dock Assembly

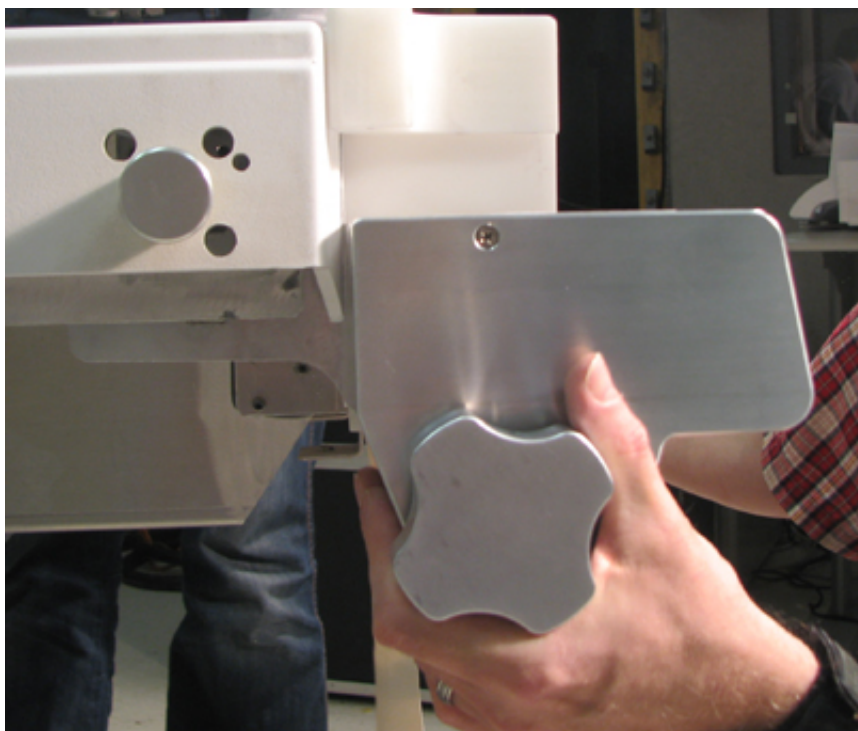

6.3.2 Replacement Dock Assembly Installation

1. Position the replacement dock assembly into the end of the table.
2. Secure the dock with the two-inch wide duct tape.
3. Reinstall the six Allen-head screws (three on each side) to mount the dock to the table. Do not tighten the screws at this time.
4. Move the cradle back to the Home position to prepare to move the table.
5. Move the OR-Compatible Table into the Operating Room.
6. Make sure the table is set to the highest position by pumping the Up pedal on the table.
7. Carefully dock the OR-Compatible Table to the Maquet Surgical Table.
8. Verify that the OR-Compatible Table is close enough to the Maquet Surgical Table to allow for undocking of the patient transport. The two retaining pins must be below the patient transport as shown below.

Illustration 7-22: Surgical Table Docking Adjustment

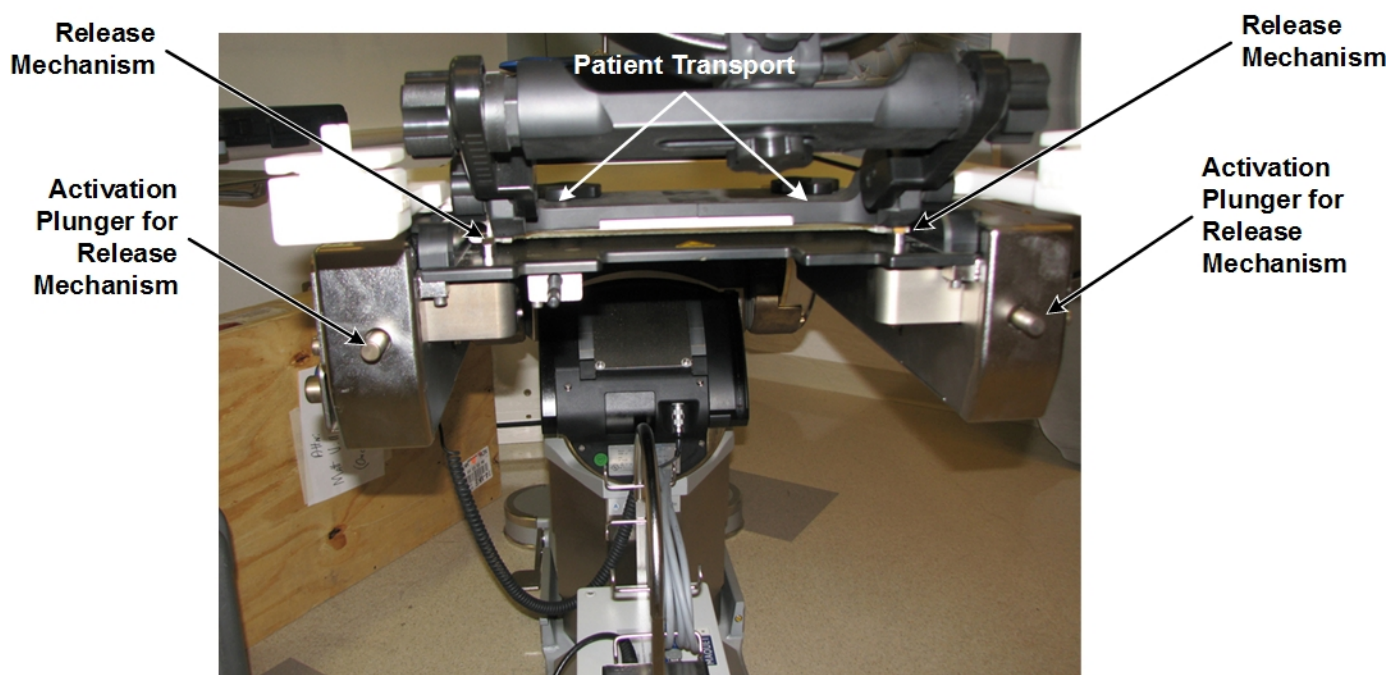

9. Ensure the two tables are parallel with each other.
10. With the tables docked, tighten the dock assembly's six Allen-head screws.
11. Undock and re-dock the tables to ensure proper operation and alignment. Make adjustments to the dock mounting position as necessary.
12. Move the OR-Compatible Table from the Operating Room to the service area.

13. Reinstall the Rear Side Bumpers.
14. Reinstall the Upper Rear Side Covers.
15. Reinstall the Rear Bottom Cover.

6.4 Finalization

Clean the table surfaces as necessary. Notify hospital personnel that the table is ready to be wiped down and prepared for patient use.

7 Touch and Go Strip and Side Bumper Replacement

7.1 Personnel Requirements

Personnel Requirements	Preliminary Reqs	Procedure	Finalization
1		75 mins	

7.2 Overview

On each side of the Patient Table along the top edge, there are Touch and Go (T and G) strips that run the length of the Table. A closer look at the T and G strips internally shows that each side is made up of three individual strips connected by jumpers. If one of the strips fails, you replace only *that* specific strip. In other words, you do not have to replace the entire side. Be certain of which strip you are replacing since each one has a separate part number. See the *Replacement Parts* section below for all part numbers.

NOTE: The Touch and Go feature is not included on the OR-Compatible Table.

This procedure also covers replacing the Table Side Bumpers. The convenience of replacing any of the three Side Bumpers per either side of the Table is that a Touch and Go strip is contained within a Side Bumper. So whether you need to replace a Touch and Go strip or a Side Bumper, you need to remove the Touch and Go strip from its respective Side Bumper anyway, thus one procedure covers both items.

NOTE: The OR-Compatible Table has unique part numbers for the Rear LH Side Bumper and the Rear RH Side Bumper.

7.3 Preliminary Requirements

7.3.1 Tools and Test Equipment

Item	Qty	Effectivity	Part#	Manufacturer
5/32 inch Allen wrench	1	-	-	-
3/8 inch open end wrench	1	-	-	-
Medium Phillips screwdriver	1	-	-	-

7.3.2 Replacement Parts


Item	Qty	Effectivity	Part#	Manufacturer
Front RH T and G Assembly	1	-	See FRU manual.	-
Arm Board T and G Assembly	1	-	See FRU manual.	-
Rear RH T and G Assembly	1	-	See FRU manual.	-
Front LH T and G Assembly	1	-	See FRU manual.	-

Item	Qty	Effectivity	Part#	Manufacturer
Arm Board T and G Assembly	1	-	See FRU manual.	-
Rear LH T and G Assembly	1	-	See FRU manual.	-
Front RH Side Bumper	1	-	See FRU manual.	-
Arm Board Bumper (use for RH or LH side)	1	-	See FRU manual.	-
Rear RH Side Bumper	1	-	See FRU manual.	-
Rear RH Side Bumper for OR-Compatible Table	1	-	See FRU manual.	-
Front LH Side Bumper	1	-	See FRU manual.	-
Rear LH Side Bumper	1	-	See FRU manual.	-
Rear LH Side Bumper for OR-Compatible Table	1	-	See FRU manual.	-

7.4 Procedure

1. Undock the Patient Table and move it out of the magnet room.

Illustration 7-23: Table Covers

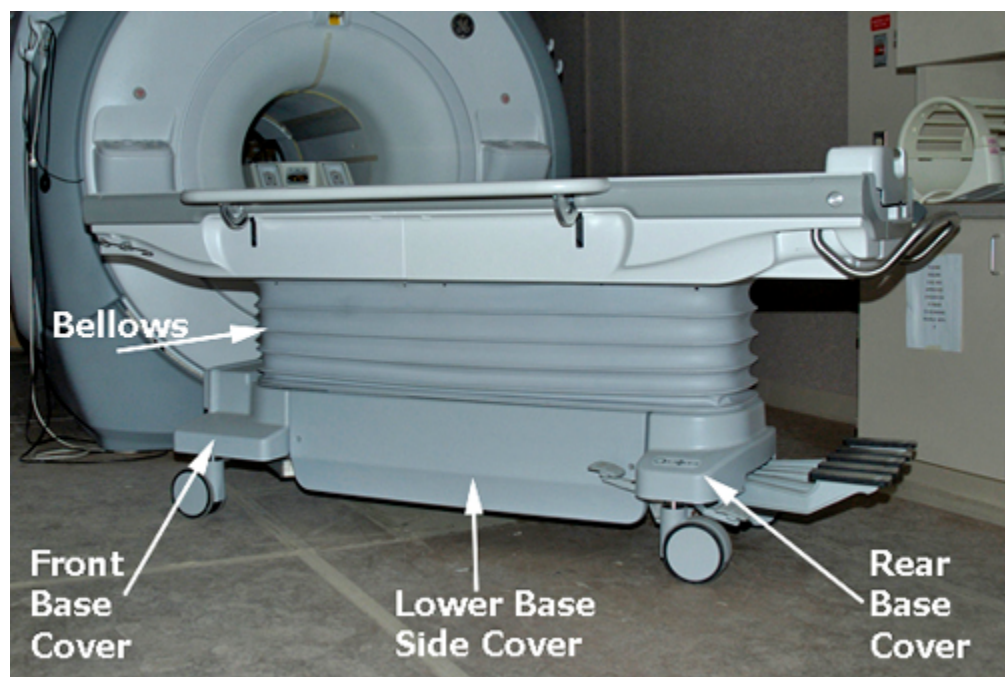


2. At the front of the Table, remove the Front Bottom Cover (4 screws). See illustration below.

Illustration 7-24: Front Bottom Cover




Illustration 7-25: Front Bottom Cover

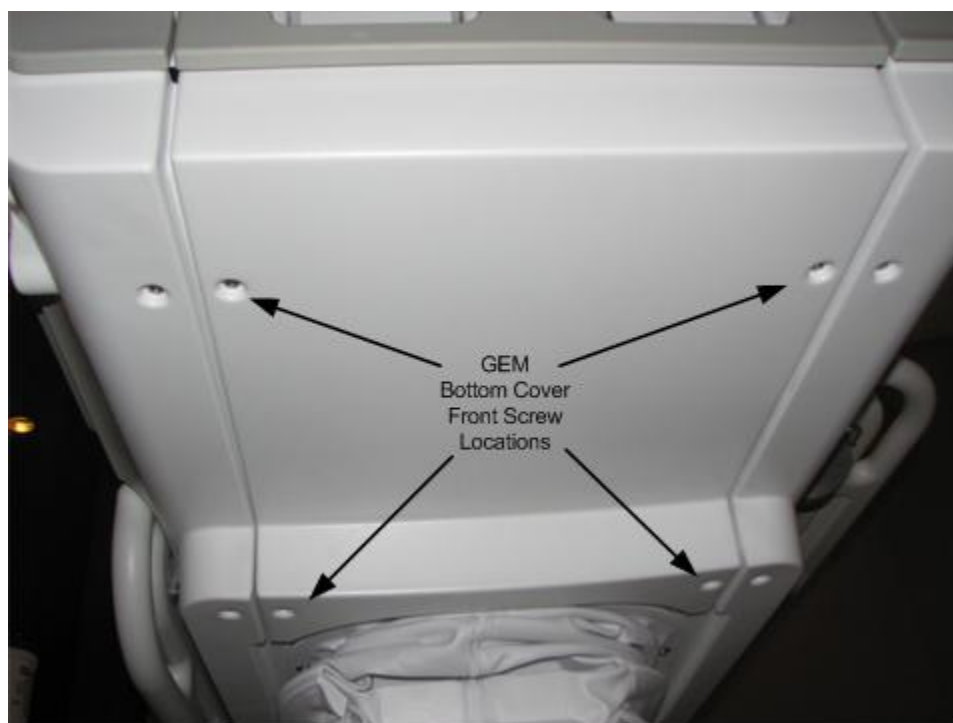

3. At the foot of the Table, remove the Rear Bottom Cover (located above the pedals, 4 screws).

4. Remove the screws from the top of the bellows on the side being worked on (five screws on each side holds the Upper Bellows to the bottom side of the Table). Do not remove the bellows from the Table; leave it in the down position.
5. (Skip this step for the OR-Compatible Table). To remove the Table handle, at the foot of the Table rotate the handle down by pushing the silver button on the left side of the Table to release the Table. Use a 5/32 inch Allen wrench to remove hex socket screws, one on each side, from the Table handle. See illustration below. Push the silver button again to rotate the handle back up to allow the handle to slip off.

Illustration 7-26: Table Handle in Down Position

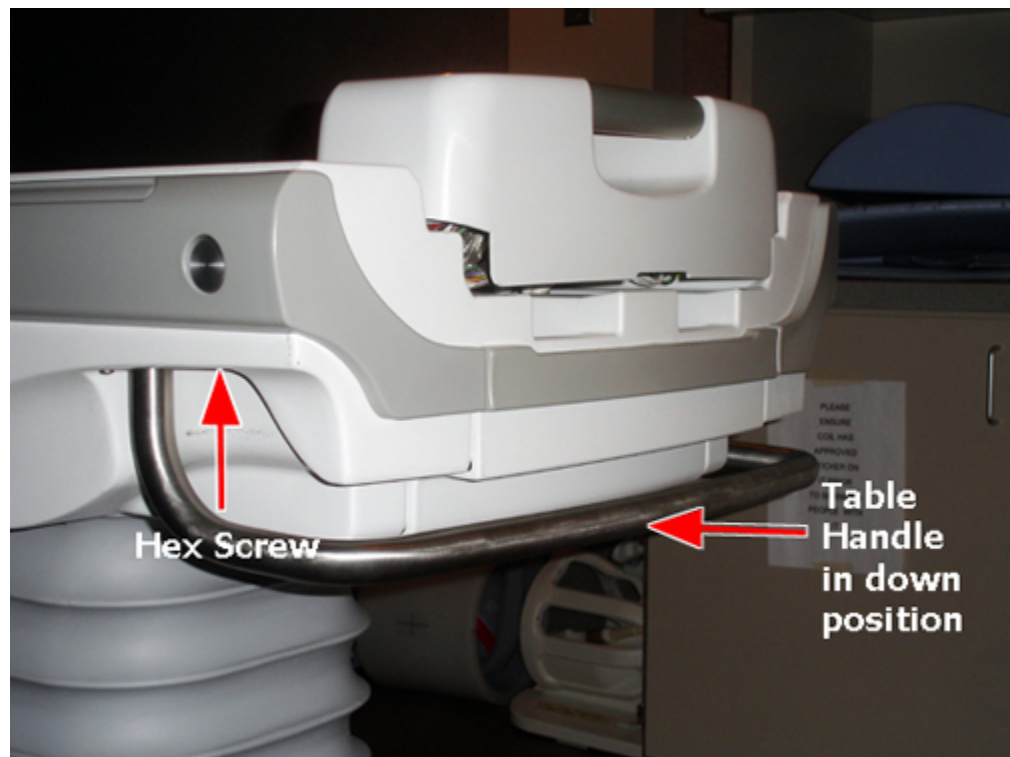


- NOTE: Side Cover Removal Strategy.** determine which T and G Strip or Side Bumper to replace, and then pay attention to the following:
- If replacing the front T and G Strip or Side Bumper, then remove the front Table Side Cover (4 screws).
 - If replacing the middle T and G Strip or Side Bumper, then remove both Front and Rear Table Side Covers (8 screws).
 - If replacing the rear T and G Strip or Side Bumper, then remove the rear Table Side Cover (4 screws).
6. On the bottom right or left side of the Table, remove appropriate Side Cover screws and then slide the Side Cover down. See illustration below.

Illustration 7-27: Side Cover, Bottom Side Screws




Illustration 7-28: Side Cover, Bottom Side Screws

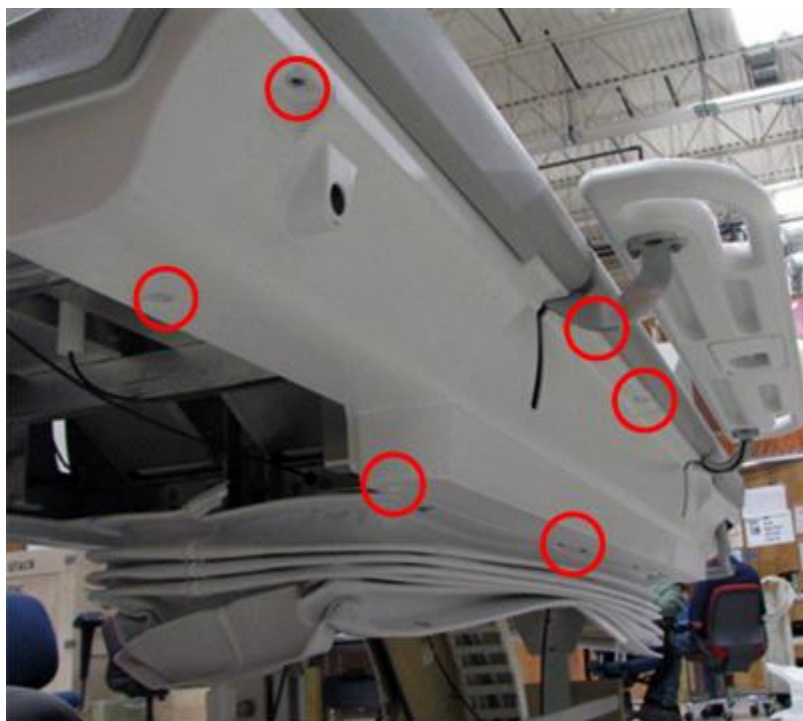

7. (Skip this step for the OR-Compatible Table.) Disconnect the connection(s) from the touch and go strip/side bumper assembly that is being replaced.
8. Use a 3/8 inch open end wrench to remove nuts and washers from the appropriate front, middle, or rear T and G Strip/Side Bumper studs.

Illustration 7-29: Side Cover, Side Bumper, T and G Strip Diagram

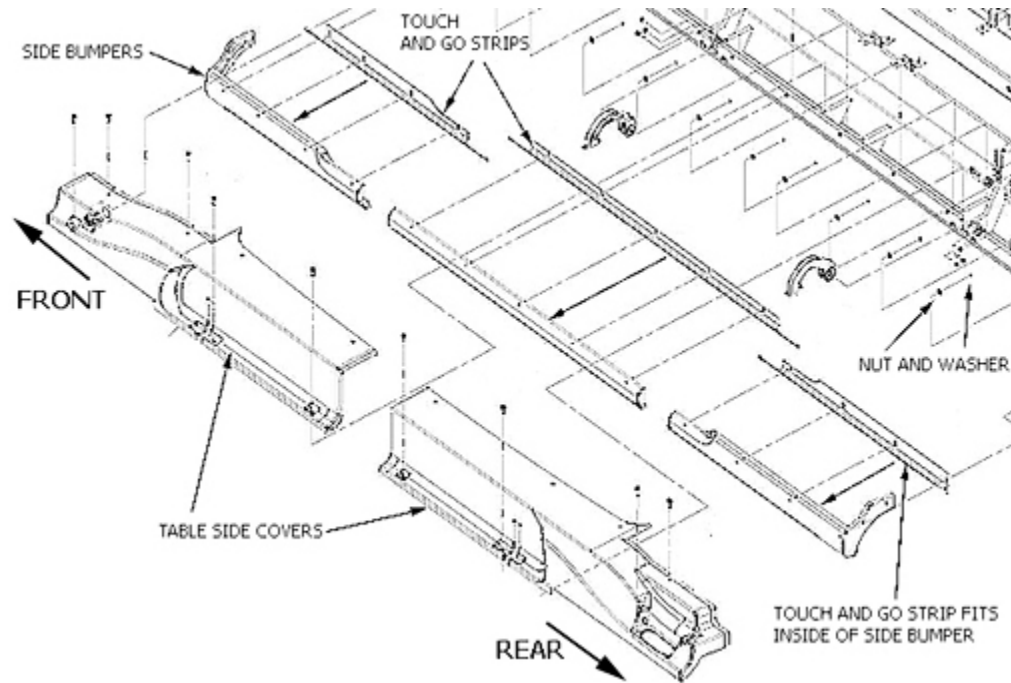

9. Slide the T and G Strip/Side Bumper assembly away from the Table chassis. See [Illustration 7-29](#).
10. Disassemble the T and G Strip from the Side Bumper and replace desired part.
11. Reverse Steps 7-10 to reassemble T and G Strip(s), Side Bumper(s), and Table Side Cover(s).
12. (Skip this step for the OR-Compatible Table.) Dock the Table and ensure that the touch and go strips are working correctly.
13. Reinstall the Table side covers.
14. Replace Table Handle screws, one on each side using a 5/32 inch Allen wrench.
15. Reattach the top side Bellows screws to bottom side of Table.

7.5 Finalization

No finalization steps.


This page left intentionally blank.

Chapter 8 Troubleshooting

1 Table Down Inhibit Not Working When Docked

The Down Pedal on the OR-Compatible Table is designed to only lower the table when the table is not docked to the scanner. See [Illustration 8-1](#) to [Illustration 8-4](#) for operation of the Table Down Inhibit mechanism, located on the right side when docked to the scanner.

Illustration 8-1: Table Docked: Down Pedal Not Depressed

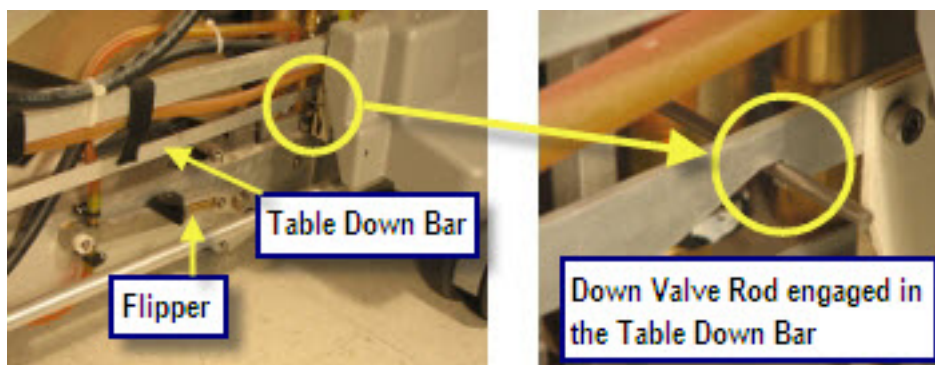

When the table is docked, depressing the Down Pedal results in the Flipper rotating clockwise to raise the Table Down Bar. See illustrations, below.

Illustration 8-2: Table Docked; Down Pedal Depressed



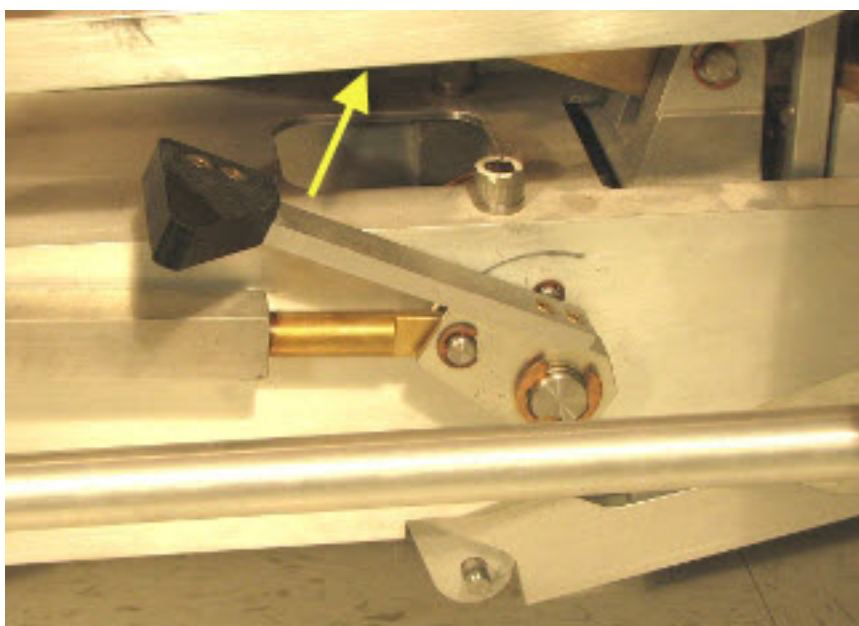


Illustration 8-3: Flipper Raising the Table Down Bar

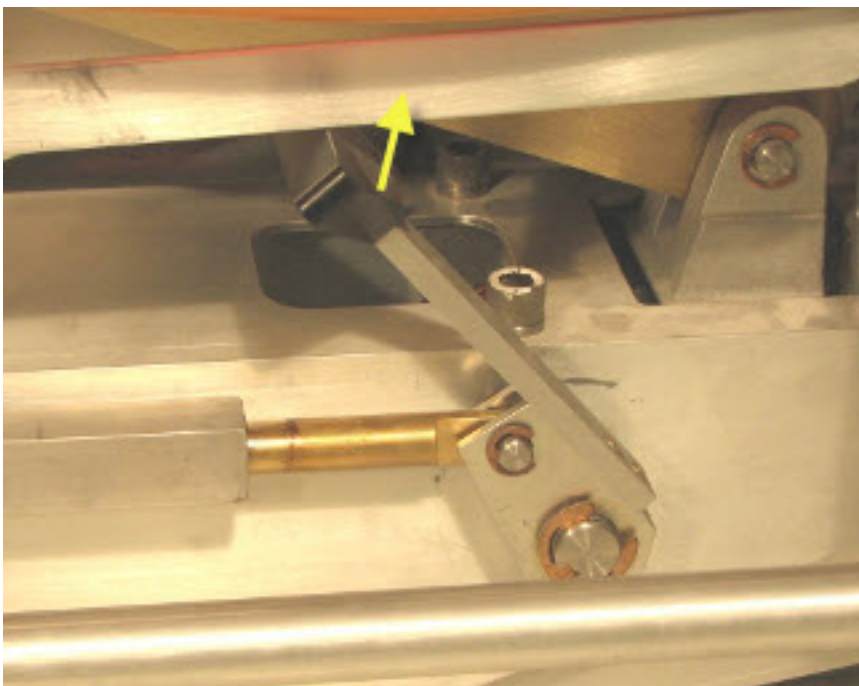

When raised by the Flipper, the Table Down Bar prevents the Down Valve Rod from being pulled forward to lower the table. See the illustration below.

Illustration 8-4: Flipper Raising Table Down Bar and Inhibiting Down Valve Rod

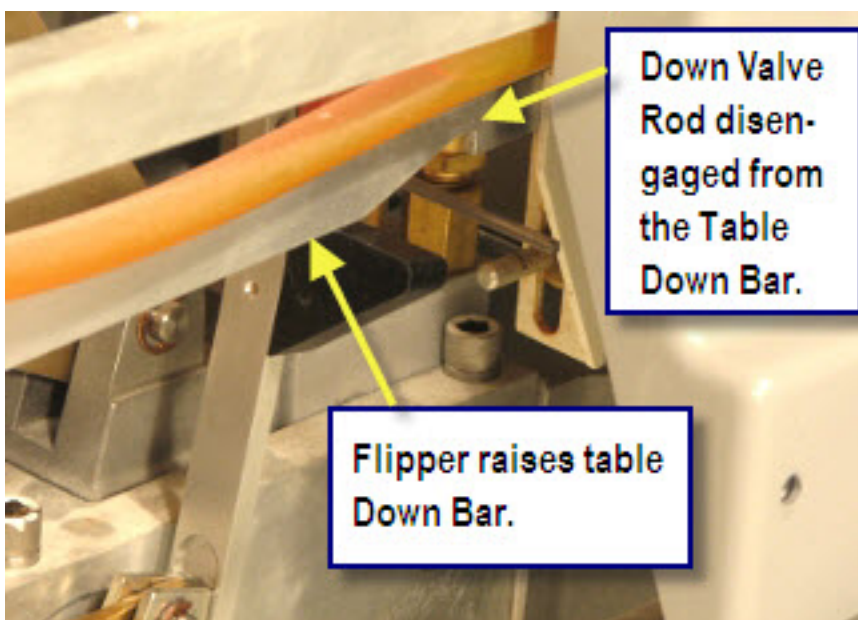

Causes

- Down Valve Rod is bent up.

- Loose screws on the Flipper Assembly.

Solution

Confirm that the Table Down Rod was not bent during other servicing activities. Refer to Cradle Home Table Down Interlock Adjustment for Table Down Rod procedures.

This page left intentionally blank.

© 2010 General Electric Company.

GE Medical Systems, a General Electric Company, going to market as GE Healthcare.

www.gehealthcare.com

