

GE Healthcare

Signa HD, HDx & HDxt 3.0T Block Diagrams & Supplemental Schematics (Release 12.x, 14.x, 15.x & 16.x)

Technical Publication



5133320
Rev 12

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation “**damage in shipment**” written on **all** copies of the freight or express bill **before** delivery is accepted or “signed for” by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

File a report with

- Call 1–800–548–3366 and use option 8.
- Contact your local service coordinator for more information on this process.

Rev. 08/15/2003

© 2004 – 2011, General Electric Company
GE Medical Systems, a General Electric Company, going to market as GE Healthcare
3200 N. Grandview Boulevard
Waukesha, Wisconsin 53188
USA

www.gehealthcare.com



WARNING

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

AVERTISSEMENT

- CE MANUEL DE MAINTENANCE N'EST DISPONIBLE QU'EN ANGLAIS.
- SI LE TECHNICIEN DU CLIENT A BESOIN DE CE MANUEL DANS UNE AUTRE LANGUE QUE L'ANGLAIS, C'EST AU CLIENT QU'IL INCOMBE DE LE FAIRE TRADUIRE.
- NE PAS TENTER D'INTERVENTION SUR LES ÉQUIPEMENTS TANT QUE LE MANUEL SERVICE N'A PAS ÉTÉ CONSULTÉ ET COMPRIS.
- LE NON-RESPECT DE CET AVERTISSEMENT PEUT ENTRAÎNER CHEZ LE TECHNICIEN, L'OPÉRATEUR OU LE PATIENT DES BLESSURES DUES À DES DANGERS ÉLECTRIQUES, MÉCANIQUES OU AUTRES.

WARNUNG

- DIESES KUNDENDIENST-HANDBUCH EXISTIERT NUR IN ENGLISCHER SPRACHE.
- FALLS EIN FREMDER KUNDENDIENST EINE ANDERE SPRACHE BENÖTIGT, IST ES AUFGABE DES KUNDEN FÜR EINE ENTSPRECHENDE ÜBERSETZUNG ZU SORGEN.
- VERSUCHEN SIE NICHT, DAS GERÄT ZU REPARIEREN, BEVOR DIESES KUNDENDIENST-HANDBUCH ZU RATE GEZOGEN UND VERSTANDEN WURDE.
- WIRD DIESE WARNUNG NICHT BEACHTET, SO KANN ES ZU VERLETZUNGEN DES KUNDENDIENSTTECHNIKERS, DES BEDIENERS ODER DES PATIENTEN DURCH ELEKTRISCHE SCHLÄGE, MECHANISCHE ODER SONSTIGE GEFAHREN KOMMEN.

AVISO

- ESTE MANUAL DE SERVICIO SÓLO EXISTE EN INGLÉS.
- SI ALGÚN PROVEEDOR DE SERVICIOS AJENO A GEMS SOLICITA UN IDIOMA QUE NO SEA EL INGLÉS, ES RESPONSABILIDAD DEL CLIENTE OFRECER UN SERVICIO DE TRADUCCIÓN.
- NO SE DEBERÁ DAR SERVICIO TÉCNICO AL EQUIPO, SIN HABER CONSULTADO Y COMPRENDIDO ESTE MANUAL DE SERVICIO.
- LA NO OBSERVANCIA DEL PRESENTE AVISO PUEDE DAR LUGAR A QUE EL PROVEEDOR DE SERVICIOS, EL OPERADOR O EL PACIENTE SUFRAN LESIONES PROVOCADAS POR CAUSAS ELÉCTRICAS, MECÁNICAS O DE OTRA NATURALEZA.

ATENÇÃO

- ESTE MANUAL DE ASSISTÊNCIA TÉCNICA SÓ SE ENCONTRA DISPONÍVEL EM INGLÊS.
- SE QUALQUER OUTRO SERVIÇO DE ASSISTÊNCIA TÉCNICA, QUE NÃO A GEMS, SOLICITAR ESTES MANUAIS NOUTRO IDIOMA, É DA RESPONSABILIDADE DO CLIENTE FORNECER OS SERVIÇOS DE TRADUÇÃO.
- NÃO TENHA TENTADO REPARAR O EQUIPAMENTO SEM TER CONSULTADO E COMPREENDIDO ESTE MANUAL DE ASSISTÊNCIA TÉCNICA.
- O NÃO CUMPRIMENTO DESTA AVISO PODE POR EM PERIGO A SEGURANÇA DO TÉCNICO, OPERADOR OU PACIENTE DEVIDO A CHOQUES ELÉTRICOS, MECÂNICOS OU OUTROS.

AVVERTENZA

- IL PRESENTE MANUALE DI MANUTENZIONE È DISPONIBILE SOLTANTO IN INGLESE.
- SE UN ADDETTO ALLA MANUTENZIONE ESTERNO ALLA GEMS RICHIEDE IL MANUALE IN UNA LINGUA DIVERSA, IL CLIENTE È TENUTO A PROVVEDERE DIRETTAMENTE ALLA TRADUZIONE.
- SI PROCEDA ALLA MANUTENZIONE DELL'APPARECCHIATURA SOLO DOPO AVER CONSULTATO IL PRESENTE MANUALE ED AVERNE COMPRESO IL CONTENUTO.
- NON TENERE CONTO DELLA PRESENTE AVVERTENZA POTREBBE FAR COMPIERE OPERAZIONI DA CUI DERIVINO LESIONI ALL'ADDETTO ALLA MANUTENZIONE, ALL'UTILIZZATORE ED AL PAZIENTE PER FOLGORAZIONE ELETTRICA, PER URTI MECCANICI OD ALTRI RISCHI.

警告

- ・ このサービスマニュアルには英語版しかありません。
- ・ GEMS以外でサービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。
- ・ このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないで下さい。
- ・ この警告に従わない場合、サービスを担当される方、操作員あるいは患者さんが、感電や機械的又はその他の危険により負傷する可能性があります。

注意：

- 本维修手册仅存有英文本。
- 非 GEMS 公司的维修员要求非英文本的维修手册时，客户需自行负责翻译。
- 未详细阅读和完全了解本手册之前，不得进行维修。
- 忽略本注意事项会对维修员，操作员或病人造成触电，机械伤害或其他伤害。

REVISION HISTORY

REV	DATE	PRIMARY REASON FOR CHANGE
1. . .	Oct 28, 2004	Initial Release
2 . . .	May 11, 2005	Overview Tab – Sec 3: Added new acronyms and descriptions. System Tab – Added XW8200 host computer, removed CAN Link for ACGD, added HFD–S CAN Link, and other minor corrections. Operator Workspace Tab – Added XW8200 host computer and other minor corrections. Patient Handling Tab – Made minor corrections. Gradient Tab – Removed ACGD cabinet pages and inserted HFD–S Cabinet page. MGD/RRF/RF Tab – Corrected minor errors and added informational callouts. Added MRC II Tool & Cable and MRC II Tool – 8 Channel Connector pages.
3 . . .	Oct 15, 2005	System Tab – Added J116 and J119 to TAC Cabinet and relocated Terminator. Patient Handling Tab – Added test points and LED's to Magnet Enclosure and PHPS. Added MNS connection information to SSM. MGD/RRF/RF Tab – Added LED info to Exciter and Receiver boards. Corrected RRF Receive Chain connection labels. Added MNS connection info to Driver Module and SSM. Added MC Preamp Bias Filter pinouts. Added LED and Pinout info for 16 Ca'hannel Switch.
4 . . .	May 10, 2006	Major update to manual with the addition of HDx (Release 14.x) product. Changed title page and banners on all pages. Also added 3T94 information throughout this manual. Overview Tab – Sec 3: Added new acronyms and descriptions. System Tab – Added HP9300 computer info, 32 Channel option info, and the new RFS cabinet power distribution using two power strips. Revised all sheet references. OW Tab – Added HP9300 computer info. MGD/RRF/RF Tab – Added HDx (14.x) info as needed which includes VRE modules: ICN's, 16/24 Port Ethernet Switch, Infiniband board and module. Added HDx A, B, and C port connection information for the HDx LPCA. Added MCR III Tool
5 . . .	Oct 20, 2006	Added Environmental Sensor to System/RFS cabinet in System. Operator Workspace, and MGD/RRF/RF Sections (MRIhc19337). MGD/RRF/RF Tab –Added Port 2 to Port 2 connection for 4 ICN configuration. Added 4 ICN configuration without Infiniband switch. Updated Multi–Coil Adapter by adding MC Bias 17–24 for 14.x systems. Updated several pages with additional MC signal information (MRIhc18215).
6 . . .	Dec 20, 2006	SYSTEM Tab – Added 32 Channel configuration without Infiniband Switch. OW Tab – Added Audio IN/OUT signals for PC Host for 14.x system. GRADIENT Tab – Added Pang block configuration info (MRIhc20732)
7 . . .	Jun 6, 2007	SYSTEM Tab – Added note for designating components required for 16 channel systems. MGD/RRF/RF Tab –Added notes to designate parts not included with newer 8 Channel systems and for parts that are required for 16 Channel systems.
8 . . .	Oct 26, 2007	SYSTEM Tab – Added CAM Chassis designation to MGD Chassis where applicable. Added CAN LINK block diagram for RFS cabinet with CAM Chassis. MGD/RF/RRF Tab – Added board slot locations for CAM Chassis.
9 . . .	Sep 11, 2008	Changed title of manual and all page headers to <i>Signa HD, HDx & HDxt 3.0T (12.x, 14.x & 15.x) Block Diagrams & Supplemental Schematics</i> . Updates for HDxT (15.x) in all sections. Added information on Sun ICN's.. Added depopulated 16–channel switch for 8–channel installations.
10 . .	Dec 17, 2008	SYSTEM Tab – Added new environmental sensor diagram (ECO 2055310). PATIENT HANDLING Tab – Revised E–Stop circuit voltage from 24V to 12V in Teal PDU (PQR 13207104).
11	Jun 30, 2010	SYSTEM Tab – Added new Z400 computer page (ECO2056587) and added 5337894–3 and –4 ICN's (ECO 2101939). MGD/RRF/RF Tab – Added new ICN's (ECO 2101939).
12	Jan 20, 2011	Update entire manual headers and title page to include 16.x software version (DOC0853690). Added XW8400 Host computer to SYSTEM & OW tabs. Removed terminator from all SCSI towers. Removed J11, SCSI tower power cord, from power strip for Z400 host computers (MRIhc54665).

LIST OF EFFECTIVE PAGES

PAGE	REV	PAGE	REV	PAGE	REV	PAGE	REV	PAGE	REV
Title Page	12	1–1	12	OPERATOR		GRADIENT		MGD/RRF/RF	
GE Logo Page	–	2–1	12	WORKSPACE		Tab MR418	–	Tab MR451	–
2128126	1*	3–1 to 3–6	12	Tab 2158195	–	1–1 to 1–3	12	i to ii	12
A	12			i	12			1–1 to 1–47	12
i	12			1–1 to 1–9	12	PDU		2–1	12
		SYSTEM		PATIENT HANDLING		Tab MR406	–	COOLING SYSTEM	
OVERVIEW		Tab MR364	–	Tab MR405	–	i	12	Tab 2308173	–
Tab MR367	–	i	12	i	12			i	12
i	12	1–1 to 1–21	12	1–1 to 1–9	12			1–1 to 1–2	12

* This revision number/letter corresponds to the indicated document's revision control system.

CONTENTS

OVERVIEW

Purpose
Standard Conventions
Abbreviations



SYSTEM

System Block Diagrams



OPERATOR WORKSPACE

Operator Workspace Block Diagrams



PATIENT HANDLING

Patient Handling Block Diagrams



GRADIENT

Gradient Block Diagrams (HFD & ACGD)



PDU

PDU Block Diagrams

Note:

Phoenix PDU found in Vendor Manuals only. (SEE Service Methods CDRom)



MGD/RRF/RF

MGD/RRF/RF Block Diagrams
RF Supplemental Schematics



COOLING SYSTEM

Water Chillers Block Diagrams

Note:

For all Chillers and Cooling System variations (SEE Vendor Manuals on MR Service Methods CDRom)



OVERVIEW

TABLE OF CONTENTS


<u>SECTION</u>	<u>PAGE</u>
1 – PURPOSE	1-1
2 – STANDARD CONVENTIONS	2-1
3 – ABBREVIATIONS	3-1

SECTION 1 – PURPOSE

The purpose of this document is to provide an overview of the hardware which makes up the MR Signa System, and a quick means to troubleshoot major signals within the system.

It is recommended that you page through this manual to become familiar with its content and organization, beginning with System Tab, OVERALL SYSTEM. This section shows how the major subsystems are grouped by cabinet, and the main communication lines between them.

The remaining block diagrams are organized into sections which correspond to major subsystems. Some diagrams show hardware from another subsystem in order to complete an entire signal flow route. For such exceptions, the board/module title of the circuitry shown will generally help identify the applicable subsystem. See Illustration 1–1 for tab organization of this manual.

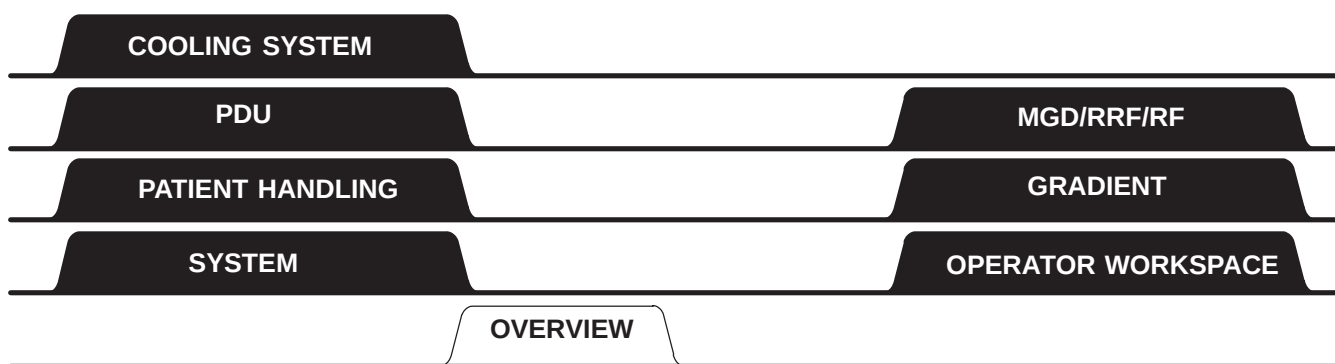


TAB ORGANIZATION FOR THIS DIRECTION
ILLUSTRATION 1

These block diagrams, along with provided diagnostics, are intended to be the primary tools for troubleshooting the MR System from a system or subsystem level on down to the specific module or board level. Since some of the RF subsystem block diagrams are presently insufficient for troubleshooting to the FRU (Field Replaceable Unit) level, schematics have been provided to supplement those block diagrams. When using the diagnostics, the block diagrams can provide a better understanding of what circuitry the diagnostics are testing, and how to further isolate a problem if required.


SECTION 2 – STANDARD CONVENTIONS

The following standard conventions are used within the block diagrams:

- The type of line used to enclose a major box on a block diagram indicates the following:

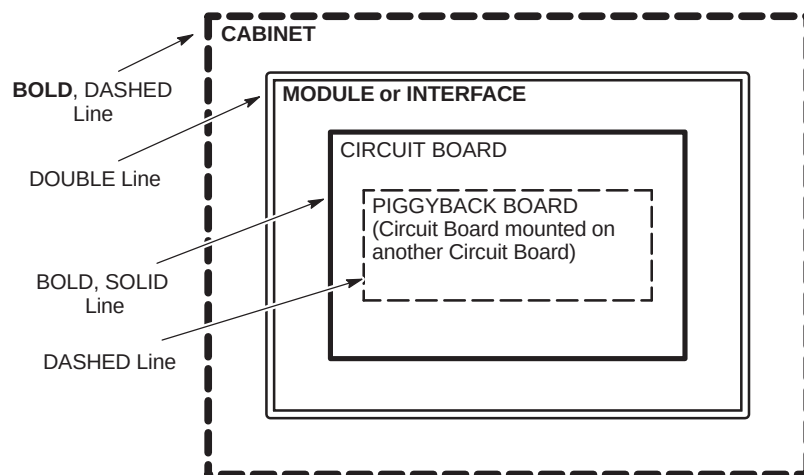

- The component designator for cabinets, modules, circuit boards, etc. is enclosed in "()" after or below their name on the block diagram.
- Logic levels: A "0" indicates a logic "low" level (typ. <0.8VDC), a "1" indicates a logic "high" level (typ. >2.0VDC).
- An "*" is used to identify an "active low" signal (e.g. "INT*" signal when at a logic "0" indicates that an interrupt is present; a logic "1" level indicates that no interrupt is present).
- A "#/" indicates the number of wires in a cable (e.g. $\frac{9}{1}$ indicates 9 wires).
- Signal flow is generally from left to right on a block diagram. Arrows are used to identify exceptions to this rule and to aid in identifying input/output signal lines on the diagram.
- The "x" symbol (e.g. —x) shown on boxes in a block diagram are used to indicate that the function within the box is enabled by the "active" signal on this input line.

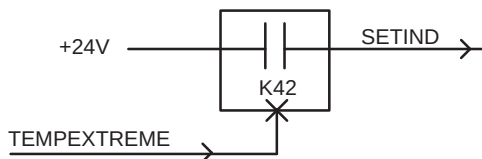

For this example, when signal TEMPEXTREME goes active ("1"), relay K42 will close (become enabled) and the signal SETIND now becomes active to set a temperature indicator.

SECTION 3 – ABBREVIATIONS

Note

This list has been compiled from a wide variety of sources. Not all of these abbreviations will apply to your system. For an electronic copy of over 1000 acronyms, see *Direction 2124201–2*.

A

AC = Alternating Current
ACGD = Advanced Concept Gradient Driver
ACK = Acknowledge
ACQ = Acquisition
A/D = Analog-to-Digital
ADC = Analog to Digital Converter
ADD or ADDR = Address
ADJ = Adjust
ADV = Advance
AGC = Automatic Gain Control
AGP = Applications Gateway Processor (Board)
AI = Applications Interface, Artificial Intelligence
ALGN = Align
ALU = Arithmetic Logic Unit
AM = Amplitude Modulation
AMP = ampere, amplifier
A/N = Alpha-numeric
ANSI = American National Standards Institute
ANT = Antenna
AP = Array Processor
APM = Analog Power Monitor
APS = Auto Pre-Scan or Acquisition Processing System (Board)
APSGT = Auto Pre-Scan Gain Control
ARCH = Archive/Remove Process
ASC = Automated Support Center or Amplifier Support Controller
ASCII = American Standard Code for Information Interchange
ASM or ASSY = Assembly
ASM = Analog Service Module
ATTEN = Attenuator
AUTO = Automatic
AUX = Auxiliary

B

B = magnetic field
BATT = Battery
BB = Buffered Data or Broadband
BC = Body Coil

BD = Board
BDS = Buffered Data Strobe
BLK = Black
BNC = BNC Connector
BRM = Body Resonance Module
BRT = Brightness
BUF = Buffer
B/W = BandWidth

C

C = centigrade, Celsius
CAB = Cabinet
CAL = Calibration
CALC = Calculate
CAN = Controller Area Network
CATTEN = Course Attenuation
CB = Circuit Breaker
CCC = CAN Core Controller card
CCW = Counter Clockwise
CD-ROM = Compact Disk Read-Only Memory
CE = Chip Enable
CF = Center Frequency
CFA = Center Frequency Adjust
CHAN or CHNL = Channel
CHK = Check
CKT = Circuit
CLR = Clear
CMD = Command
CMMR = Common Mode Rejection Ratio
CMOS = Complimentary Symmetry Metal Oxide Semiconductor
CMP = Compare or Compressor
CNT = Contrast
COAX = Coaxial
COF = Cut Off Frequency
CONFIG = Configuration
CONT = Continuous
CONV = Converter
CPD = Communication Pin Driver
CPU = Central Processing Unit

CRT = Cathode-Ray Tube
 CTL = Clock Time Control
 CTRL = Control
 CV = Configuration Variable
 CW = Clockwise or Controller Wave
 CYL = Cylinder

D

D/A = Digital-to-Analog
 DAC = Digital to Analog Converter
 dB = Decibel
 dBm = Decibels above or below one Milliwatt
 DC = Direct Current
 DCB = Driver Control Board
 DCKSTP = Dock Stop
 DD = Dynamic Disable
 DDS = Direct Digital Synthesizer
 DDSB = Dynamic Disable Switch Box
 DEC = Decimal
 DEG = Degree(s)
 DEMUX = Demultiplexer
 DEV = device
 DIA = Diagnostic Imaging Accessories
 DIAG = Diagnostic
 DIFF = Differential
 DIG = Digit
 Dig B0 = Digital B0 Adjustment
 DIP = Dual In-line Package
 DIR = Direction
 DIS = Disable
 DM = Driver Module
 DMA = Direct Memory Access
 DMAC = Direct Memory Access Controller
 DMM = Digital Multi-Meter
 DPR or DPRAM = Dual Port Random Access Memory
 DQA = Daily Quality Assurance
 DRAM = Dynamic Random Access Memory
 DRC = Dynamic Ram Controller
 DRF = Digital Receive Filter (Board)
 DRVR = Driver
 DSP = Digital Signal Processor
 DVM = Digital Volt Meter
 DWG = Drawing

E

ECC = Extended Cursor Control or Error Checking and Correction
 ECG = Electro-Cardiogram
 EEPROM = Electrically Erasable Programmable Read Only Memory
 EFB = Envelop Feedback

EPI = Echo Planer Imaging
 e.g. = for example
 EM or EMERG = Emergency
 EMI = Electro-Magnetic Interference
 EN or ENA = Enable
 ENV = Envelope
 EOS = End Of Sequence
 EOT = End Of Travel
 EOW = End Of Waveform
 EPO = Emergency Power Off
 EPROM = Erasable Programmable Read Only Memory
 ERU = Emergency Rundown Unit
 EXCITE = Expanding Applications with Multi-channel Imaging Technology
 EXER = Exerciser

F

f = farad
 F = Fahrenheit, Fuse (designator)
 FBD = Functional Block Diagram
 FBK = Feedback
 FC = footcandle
 Fc = Center Frequency
 FDBK = Feedback
 FDO = Future Delivery Order
 FE = Field Engineer
 FET = Field-Effect Transistor
 FF or FIFO = First In – First Out
 FFD = Full Field Distortion
 FFT = Fast Fourier Transform
 FID = Free Induction Decay
 FIDAT = FIFO In Data
 FIFO or FF = First In – First Out
 FILT = Filter
 FM = Frequency Modulated
 FMI = Field Modification Instruction
 FOV = Field Of View
 FREQ = Frequency
 FRNT = Front
 FRU = Field Replaceable Unit
 ft = foot
 FT = Feet
 FTP = File Transfer Protocol

G

G = gauss
 G/A = Gradient Amplifier
 GEN = Generator
 GND = Ground
 GOC AA = Global Operator Console Audio Assembly
 GRAD = Gradient
 GRN = Green
 GRY = Grey
 GP = Gradient Processor
 GSW = Gradient Switch

H

H = henry
 H = Hexadecimal
 HART = Hardware Tracker
 HC = Head Coil
 HD = Head or High Definition
 He = Helium
 He-ALARM = Helium Alarm
 HEX = Hexadecimal
 He-WARNING = Helium Warning
 HFA = High Fidelity Amplifier
 HFD = High Fidelity Driver
 HFD-S = High Fidelity Driver – Special
 HORIZ or HORZ = Horizontal
 HPC = Hardware and Protocol Control
 HPFS = Host Platform Software
 hr = hour
 HSN = High Speed Network
 HSS (D) = High Speed Serial Data Link
 HTR = Heater
 HV = High Voltage
 HV PSU = High Voltage Power Supply Unit
 HVAC = Heating, Ventilation, & Air Conditioning
 Hz = hertz

I

I = Current or Phase or Real
 IAMP = Amplifier Current (Gradient)
 ICD = Installation Certification Document
 ICL = Coil Overcurrent (Gradient)
 ICN = Image Compute Node
 ID = inner diameter
 ID = Identity, Identifier
 i.e. = that is to say
 IEC = International Electrical Code
 IEEE = Institute of Electrical & Electronics Engineers
 I/F = Interface
 in. = inch
 INFO = Information

INS, INST, or INSTR = Instruction
 INT or INTR = Interrupt
 I/O = Input/Output
 IP = Image Processor or Plate Current (RF amp) or
 Instruction Processor or Internet Protocol
 IPA = Intermediate Power Amplifier
 IPC = Interprocess Communication
 IQ = Image Quality
 IR = Image Reject
 IRF = Interface and Remote Functions (Board)
 IRF I/O = Interface and Remote Functions
 Input/Output (Board)
 IRQ = Interrupt Request
 ISA = Interrupt Start Address
 ISI = Inter-Sequence Interrupt
 ISO = Isolator or Isolation
 IT-MGD = Image Transfer-MDG (Board)
 IT-Host = Image Transfer-Host (Board)
 IVI = Interactive Vascular Imaging

J

J = Jack Connector (designator)
 JP = Jumper

K

K = Kelvin (unit of temperature) or Kilo (Thousand)
 KB = one thousand bytes of memory
 KCal = kiloCalorie
 kg = kilogram
 kHz = kilohertz
 KVA = kilo Volt Amperes
 kW = Kilowatt

L

L or LFT = Left
 LAN = Local Area Network
 lb = pound
 LBK = Loopback
 LCA = Logic Cell Arrays
 LED = Light Emitting Diode
 LFC = Load From Cold
 LHe = Liquid Helium
 LIN = Linear
 LM or LMK = Landmark
 LMP = Lamp
 LNK = Link
 LPCA = Low Profile Carriage Assembly
 LO = Local Oscillator
 LOEP = List of Effective Pages
 LONG = Longitudinal
 LVLE = Low Voltage, Low Energy

M

m = meter
 M = mega – one million
 mA = millamps
 MAG = Magnet
 MB = megabyte
 MC = Multi-Coil
 MCD = Multi-Coil switch Driver
 MCR = Multicoil Receive (Tool or Diagnostic)
 MDS = Multidrop Serial interface
 MEM = Memory
 Mfg. = Manufacturing
 MFU = Multi-Format Unit
 MGD = Multi-Generational Data acquisition (Chassis)
 MGMNT = Management
 MHz = megahertz
 mm = millimeter
 MOD = Modifier
 MON = Monitor
 MOSFET = Metal Oxide Semiconductor Field Effect Transistor
 MOV = Metal Oxide Varistor
 MPS = Manual Pre-Scan
 MR = Magnetic Resonance
 MSG = Message
 MTR = Motor
 MUX = Multiplexer
 mV = millivolts
 MW = megawatt
 mW = milliwatt

N

N = Negative
 N₂ = Nitrogen
 n.c. = normally closed
 NEC = National Electrical Code
 Ni-Cad = Nickel Cadmium
 n.o. = normally open
 NO. = Number
 NOPROC = No Processing
 NOREC = No Reconstruction
 nS = Nanosecond

O

O₂ = oxygen
 OC = Operator's Console
 OD = outside diameter
 OM = Oxygen Monitor
 OPI = Oblique Plane Imaging
 OPTO = Optical
 ORG = Orange

OT = Overtemp
 OV = Over voltage
 OVR = Over
 OVRLD or OL = Overload
 oz = ounce
 OW = Operator Workspace

P

P = Positive
 PA = Patient or Power Amplifier
 PAC = Physiological Acquisition Controller
 PAL = Programmable Array Logic or Patient Alignment Lights
 PAN = Primary Archive Node
 PAT = Patient
 PCI = Programmable Communications Interface or Peripheral Component Interconnect
 PCT = Pulse Controller Task
 PDI = Product Delivery Instruction
 PDU = Power Distribution Unit
 PEN = Penetration
 PHPS = Patient Handling Power Supply
 PIN = Positive-Intrinsic-Negative
 PIO = Programmed Input/Output
 PLL = Phase Lock Loop
 PIXEL = Picture Element
 PM = Phase Modulation, or Planned Maintenance
 PMC = PCI Mezzanine Connector
 P/N = Positive/Negative
 PNL = Panel
 POS = Position
 POT = Potentiometer
 PP = Penetration Panel
 PPC = Power PC
 P-P = Peak to Peak
 ppm = parts per million
 PREAMP = Preamplifier
 PROC = Process, Processor
 PROG = Program, -ed, -able
 PROM = Programmable Read Only Memory or Program
 PS = Power Supply
 PSD = Pulse Sequence Database
 PSG = Pulse Sequence Generator
 PSU = Power Supply Unit
 PSS = Power Supply System Board
 PT = Patient Transport or Patient Table

PUR = Purple
 PVC = Polyvinylchloride
 PW = Pulse Width
 PWA = Printed Wiring Assembly
 PWB = printed wiring board
 PWR = Power

Q

Q = Quadrature, Quality Factor or Imaginary
 Qty = Quantity
 QUAD = Quadrature or Quadrant

R

R or RT = Right
 RAM = Random Access Memory
 RCV = Receive
 RD = Read
 RDY = Ready
 RECON = Reconstruction
 REF = Reference or Reflected
 REG = Register or Regulator
 RET or RTN = Return
 REV = Revision or Reverse
 RF = Radio Frequency
 RFAMP = Radio Frequency Amplifier
 RFA = RF Amplifier
 RFI = Radio Frequency Interference or RF Interface
 (Module)
 RFS = RF System Cabinet
 RGB = Red Green Blue
 RLA or RLY = Relay
 rms = root mean square
 ROI = Region of Interest
 ROM = Read Only Memory
 RPC = Remote Control Panel
 RRF = Remote RF (Chassis)
 RRF-DIF = Remote RF-Digital Interface board
 RSR = Register
 RST = Reset
 R/W = Read/Write
 RX = Receive
 RXD = Receive Data

S

S = Switch (designator)
 SACK = Slave Acknowledge
 SAG = Sagittal
 SAR = Specific Absorption Rate
 S/C = Superconducting
 SCA = Scan Control Assembly
 SCP = Scan Control Processor (Board)
 SCR = Silicon Controlled Rectifier

SCSI = Small Computer System Interface
 SEC = second
 SEG = Segment
 SEL = Select
 SEPI = Spin Echo Planer Imaging
 SGA = Switchable Gradient Amplifier
 SGI = Silicon Graphics Inc.
 SIG = Signal
 SII = "S-2" (type of magnet)
 SIII = "S-3" (type of magnet)
 SIMMS = Single In-line Memory Module
 SIP = Single In-line Package
 SL = Slice
 SMC = Status Monitoring and Control
 SMPTE = Society of Motion Picture and Television
 Engineers
 SN = Service Note
 SNR = Signal-to-Noise Ratio
 SPI = Serial Peripheral Interface
 SPT = System Performance Test
 SPU = Signal Processing Unit
 SRAM = Static Random Access Memory
 SRF = Solid State RF Amplifier SGA or Scalable RF
 (Cabinet)
 SRF/TRF = Sequencer Related Function / Trigger and
 Rotation Function (Board)
 SRI = Scan Room Interface
 SRV = Service
 SSM = System Support Module
 SSP = Sequence Synchronous Protocol
 SST = Small Sample Test
 STAB = Stability
 STAT = Status
 STD = Standard
 STD DEV = Standard Deviation
 STIF Bd = SRF/TRF InterFace Board
 SUB-D = Subminiature D (connector)
 SW = Switch or Software
 S/W = Software
 SYNC = Synchronous
 SYS = System
 SYSID = System Identification

T

T = tesla
TAC = Twin Accessory Cabinet
TBD = To Be Determined
TBL = Table
TEMP = Temperature
TFU = Thyristor Firing Unit
TG = Transmit Gain
TGWC = TwinSpeed Gradient Water Cooler
TNS = Transient Noise Suppressor
TP = Test Point (designator)
TR or T/R = Transmit–Receive
TRMAP = Transmit/Receive Map
TRM = Twin Coil Resonator Module
TS = Terminal Strip
TSCC = TwinSpeed Cooling Cabinet
TST = Test
TTL = Transistor to Transistor Logic
TX = Transmit
TXD = Transmit Data

U

UART = Universal Asynchronous Receiver–Transmitter
uf = microfarad
UFI = Ultra Fast Imaging
UNBLK = Unblank
uP = Microprocessor
UPLMT = Up Limit
UPM = Universal Power Monitor
USART = Universal Synchronous/Asynchronous
Receiver–Transmitter
uS or usec = microsecond
UTNS = Universal Transient Noise Suppressor (Board)
UV = Under voltage

V

VAC = Volts Alternating Current
VAR = Variable
VBUS = Bus Voltage
VCR = Video Cassette Recorder
VCNTRL = Control Voltage for 8645 Techrons
VDC = Volts Direct Current
VERT = Vertical
Vp–p = Volts Peak–to–Peak
VRE = Volume Recon Engine
VRMS = Volts, Root Mean Square
VTR = Video Tape Recorder

W

WARP = Waveform And Rotation Processor
WC = Water Chiller
WR = Write
WHT = White
WIM = Workspace Interface Module
WND = Window
WORM = Write Once Read Many
WRT = Write

X

XCVR = Transceiver
XFMR = Transformer
XMT or XMIT = Transmit


Y

YEL = Yellow

SYSTEM

TABLE OF CONTENTS

<u>SECTION</u>	<u>PAGE</u>
1 – SYSTEM	
G3 & 3T94 HD OVERALL SYSTEM – OPERATOR WORKSPACE CABINET w/8000 & TABLE (Forward Production or Pre–EXCITE to HD Upgrades)	1–1
G3 & 3T94 HD OVERALL SYSTEM – OPERATOR WORKSPACE CABINET/TABLE (EXCITE to HD Upgrades)	1–2
G3 & 3T94 HD OVERALL SYSTEM – OPERATOR WORKSPACE CABINET/TABLE (Forward Production or Pre–EXCITE to HDx & HDxt Upgrades)	1–3
G3 & 3T94 HDx & HDxt OVERALL SYSTEM – OPERATOR WORKSPACE CABINETw/9300 & TABLE (Forward Production, EXCITE, or Pre–EXCITE Upgrades)	1–4
G3 & 3T94 HDx & HDxt OVERALL SYSTEM – OPERATOR WORKSPACE CABINETw/Z400 & TABLE (Forward Production, EXCITE, or Pre–EXCITE Upgrades)	1–5
G3 HD OVERALL SYSTEM – 8 Channel	1–7
G3 HD OVERALL SYSTEM – 16 Channel	1–8
G3 HDx & HDxt OVERALL SYSTEM – 8 Channel	1–9
G3 HDx & HDxt OVERALL SYSTEM – 16 Channel	1–10
G3 HDx & HDxt OVERALL SYSTEM – 32 Channel with Infiniband Switch	1–11
G3 HDx & HDxt OVERALL SYSTEM – 32 Channel without Infiniband Switch	1–12
3T94 HD UPGRADE TO HDx & HDxt OVERALL SYSTEM–16 Channel with 2 or 4 ICN's ..	1–13
3T94 Pre–EXCITE U/G TO HDx & HDxt OVERALL SYSTEM–16 Channel with 2 or 4 ICN's	1–14
G3 & 3T94 HDx & HDxt OVERALL SYSTEM (Forward Production & Upgrades)	1–15
G3 OVERALL SYSTEM – 8 or 16 Channel (HD Forward Production or HDx & HDxt Upgrades)	1–16
G3 OVERALL SYSTEM – 8 or 16 Channel (HDx & HDxt Forward Production)	1–17
3T94 OVERALL SYSTEM – 8 or 16 Channel (HDx & HDxt Upgrades)	1–18
G3 HDx & HDxt OVERALL SYSTEM – 32 Channel	1–19
G3 & 3T94 HD and HDx & HDxt OVERALL SYSTEM – CAN Link Block Diagrams	1–20
G3 HD and HDx & HDxt OVERALL SYSTEM – 8KW MNS Option	1–21

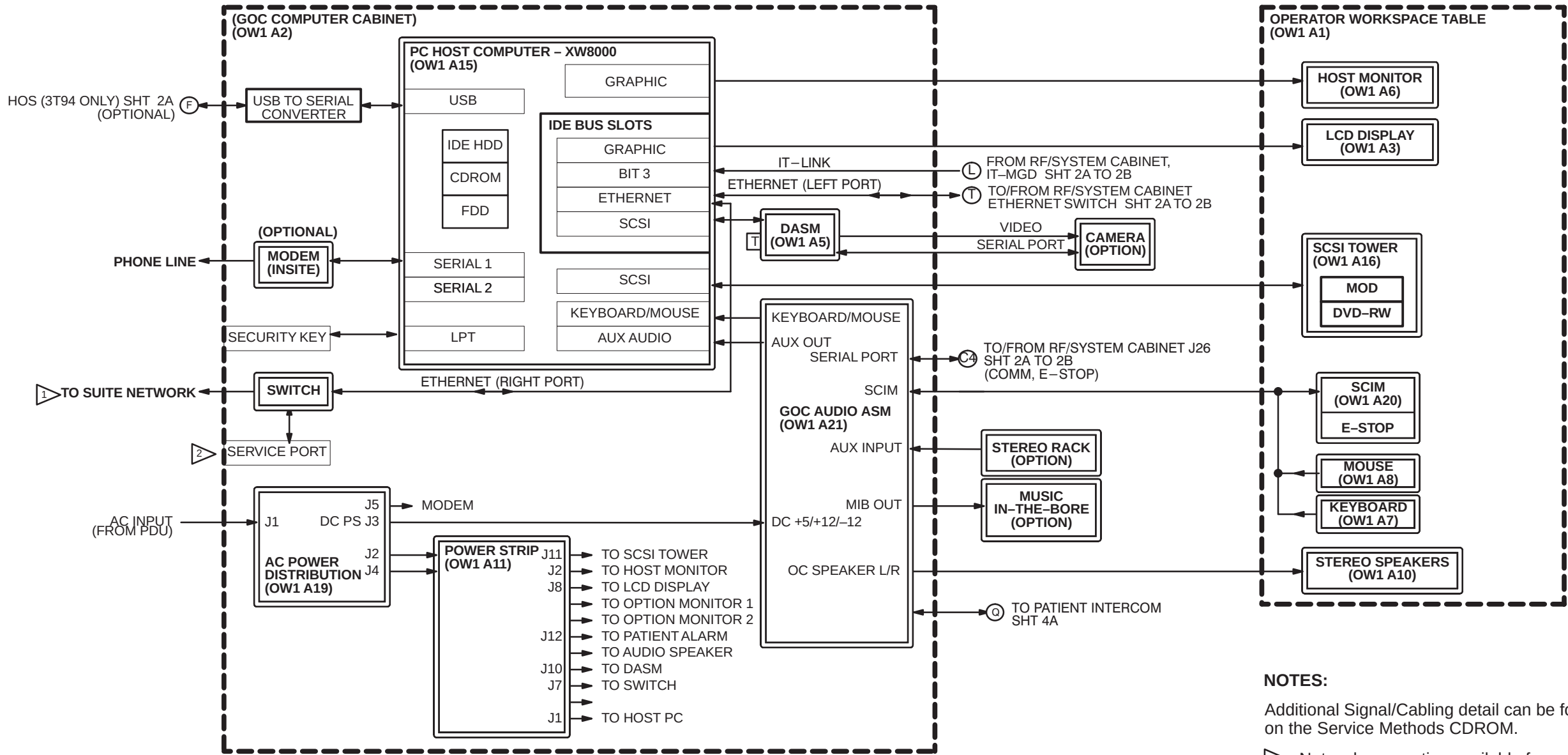

NOTES:

Additional Signal/Cabling detail can be found under Interconnects on the Service Methods CDROM.

➤ Network connection available for external suite connection


➤ Network connection available at front of goc assembly for the service laptop

Refer to the following functional block diagram for additional information:

NOTE:

1) PATIENT HANDLING

2) MGD/RRF/RF

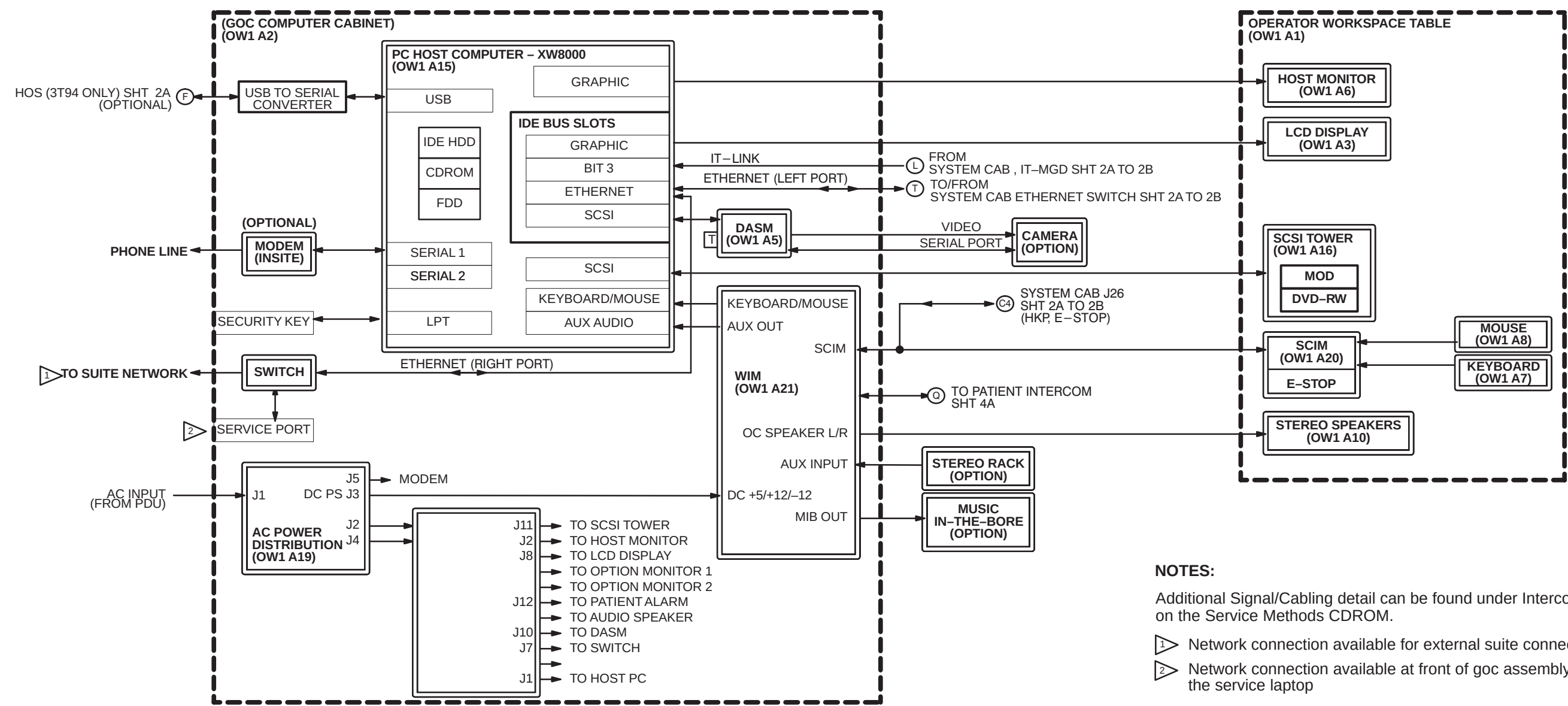


NOTES:
Additional Signal/Cabling detail can be found under Interconnects on the Service Methods CDROM.

- 1) Network connection available for external suite connection
- 2) Network connection available at front of goc assembly for the service laptop


Refer to the following functional block diagram for additional information:

NOTE:
1) PATIENT HANDLING
2) MGD/RRF/RF

G3 & 3T94 EXCITE HD (12.x) OVERALL SYSTEM
EXCITE to HD Upgrades
(OPERATOR WORKSPACE CABINET/TABLE)

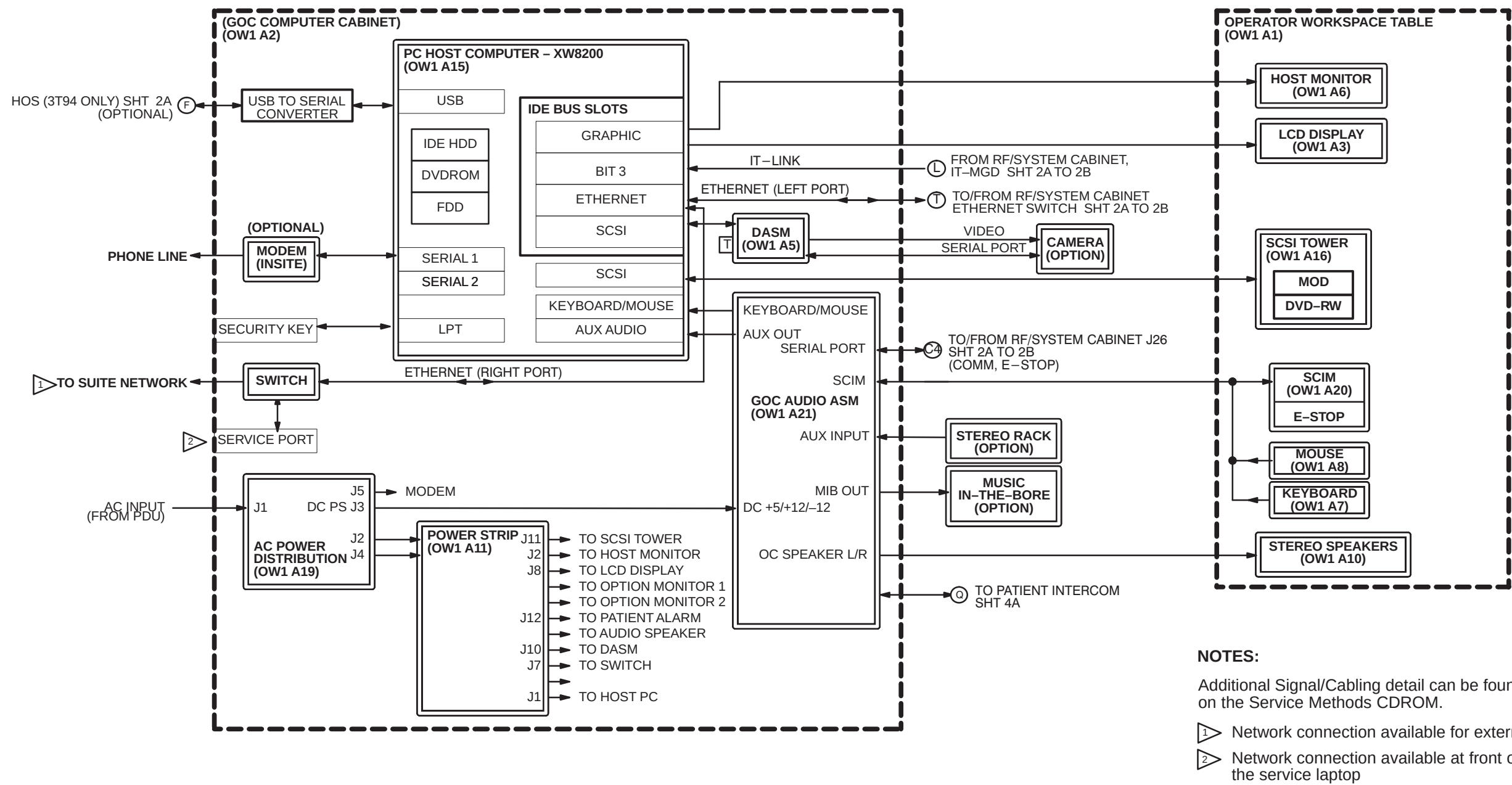


G3 & 3T94 EXCITE HD (12.x) OVERALL SYSTEM
Forward Production or Pre-EXCITE Upgrades
(OPERATOR WORKSPACE CABINET/TABLE)
SHEET 1C OF 6

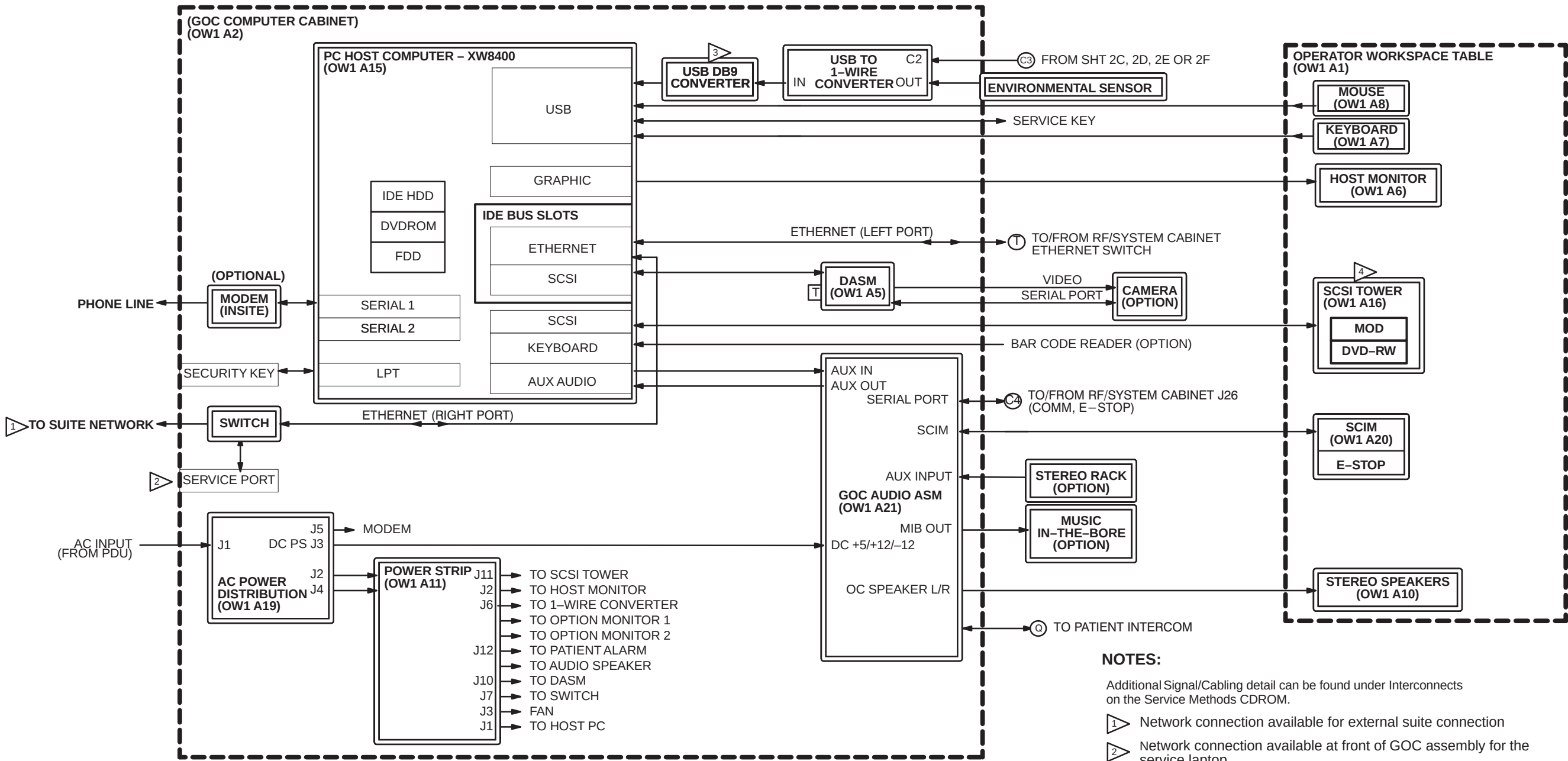


Refer to the following functional block diagram for additional information:

NOTE:
1) PATIENT HANDLING
2) MGD/RRF/RF

- NOTES:**
- Additional Signal/Cabling detail can be found under Interconnects on the Service Methods CDROM.
- 1) Network connection available for external suite connection
 - 2) Network connection available at front of GOC assembly for the service laptop
 - 3) Beginning August 2008, design changed to 1-wire to USB converter. Environmental Sensors Block Diagram utilizing 1-wire 5147136-2
 - 4) 8400 Host Computer is also available without SCSI Tower.

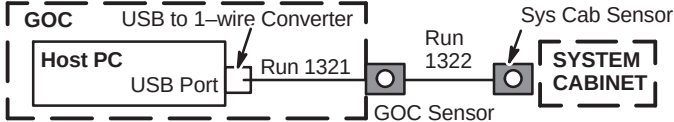

HDx (14.x) & HDxt (15.x) OVERALL SYSTEM
Forward Production, EXCITE, or Pre-EXCITE to HDx Upgrades
(OPERATOR WORKSPACE CABINET w/8400 & TABLE)

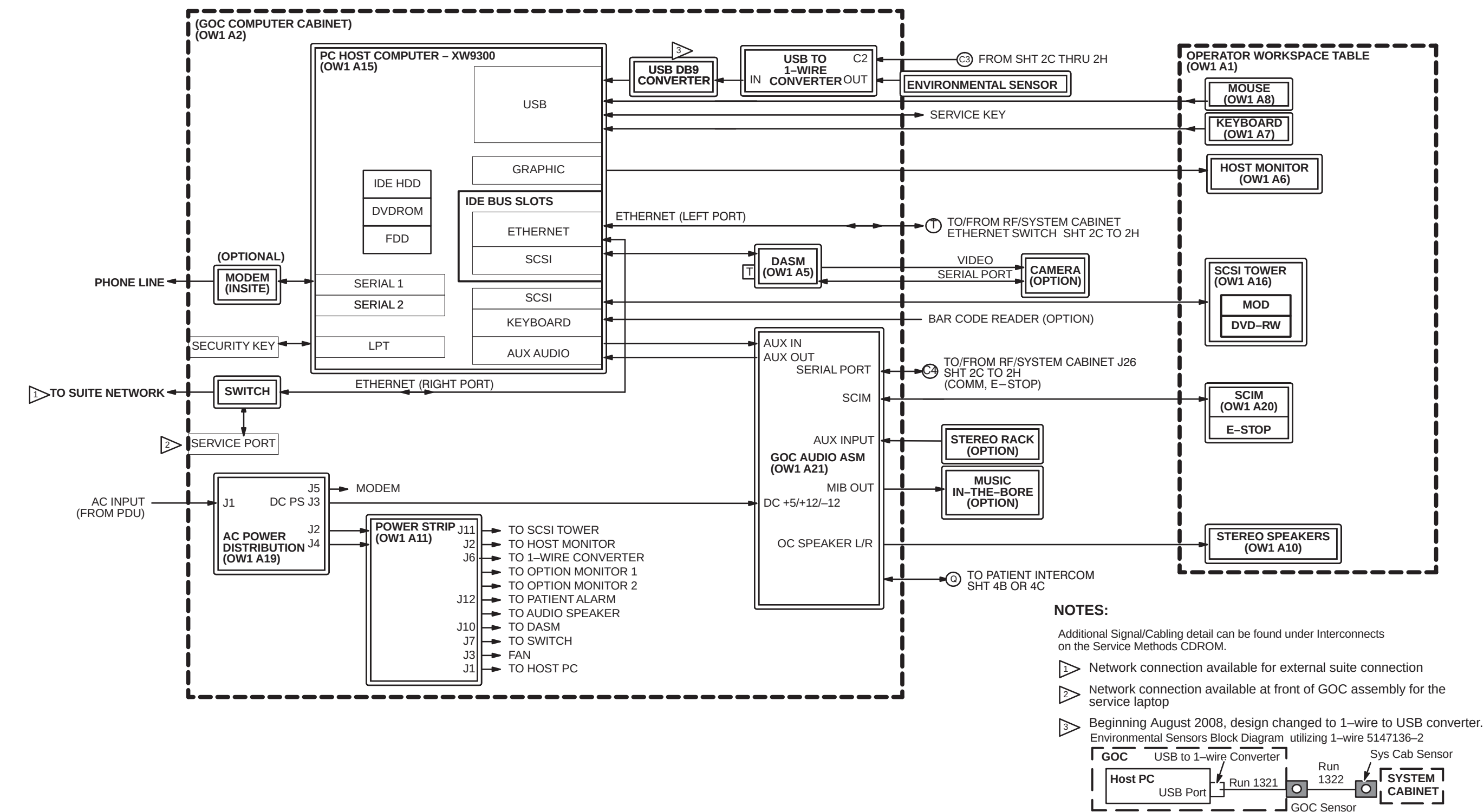


G3 & 3T94 HDx (14.x) & HDxt (15.x) OVERALL SYSTEM
Forward Production, EXCITE, or Pre-EXCITE to HDx Upgrades
(OPERATOR WORKSPACE CABINET w/9300 & TABLE)

SHEET 1E OF 6

Refer to the following functional block diagram for additional information:


NOTE:
1) PATIENT HANDLING
2) MGD/RRF/RF

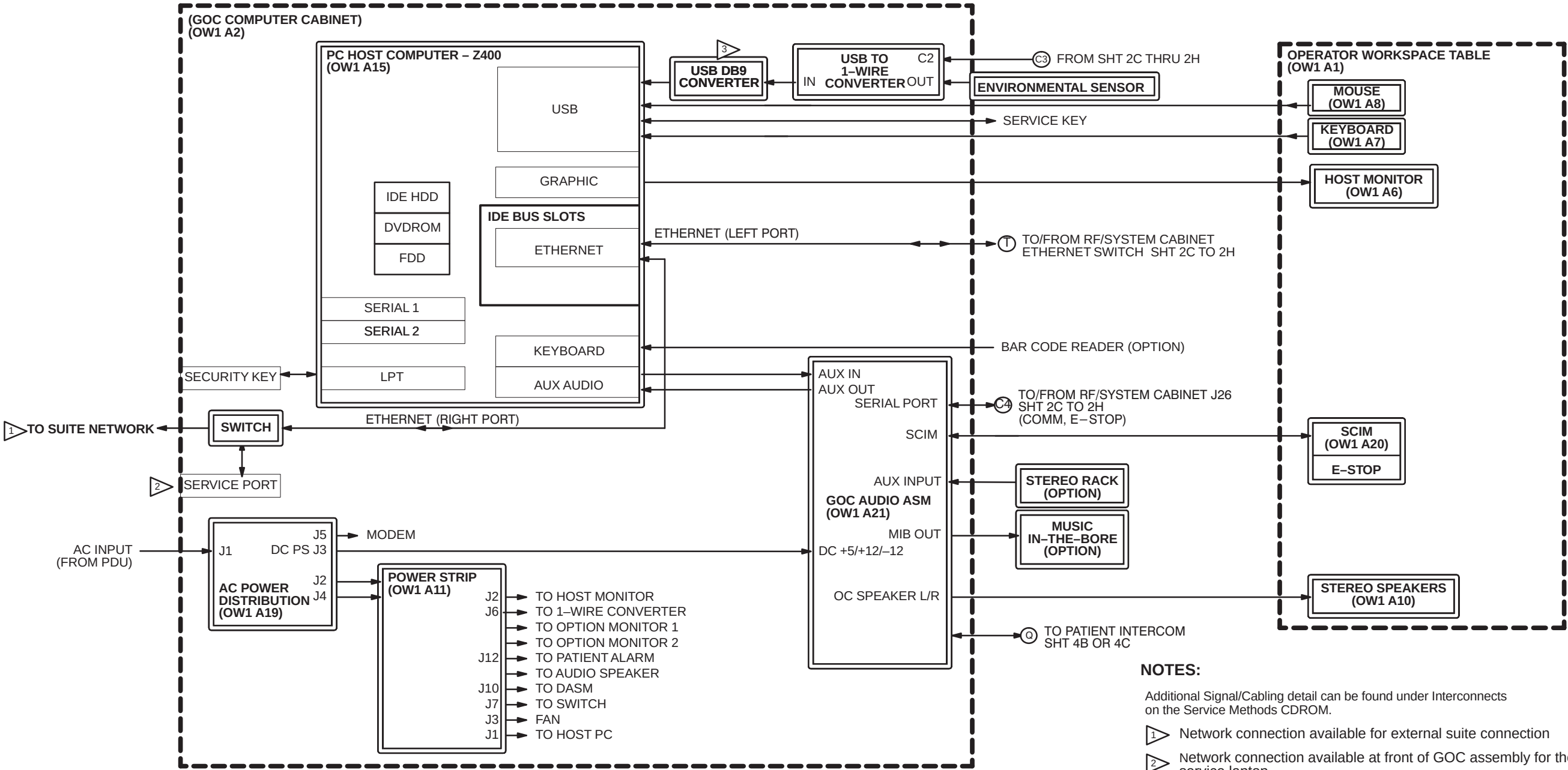

NOTES:

Additional Signal/Cabling detail can be found under Interconnects on the Service Methods CDROM.


- 1 Network connection available for external suite connection
- 2 Network connection available at front of GOC assembly for the service laptop
- 3 Beginning August 2008, design changed to 1-wire to USB converter. Environmental Sensors Block Diagram utilizing 1-wire 5147136-2

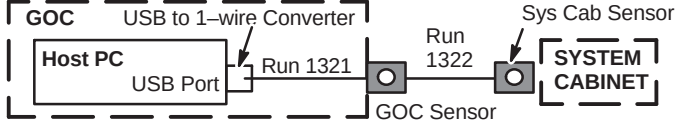


Refer to the following functional block diagram for additional information:


NOTE:
1) PATIENT HANDLING
2) MGD/RRF/RF

G3 & 3T94 HDx (14.x) & HDxt (15.x) OVERALL SYSTEM
Forward Production, EXCITE, or Pre-EXCITE to HDx Upgrades
(OPERATOR WORKSPACE CABINET w/Z400 & TABLE)

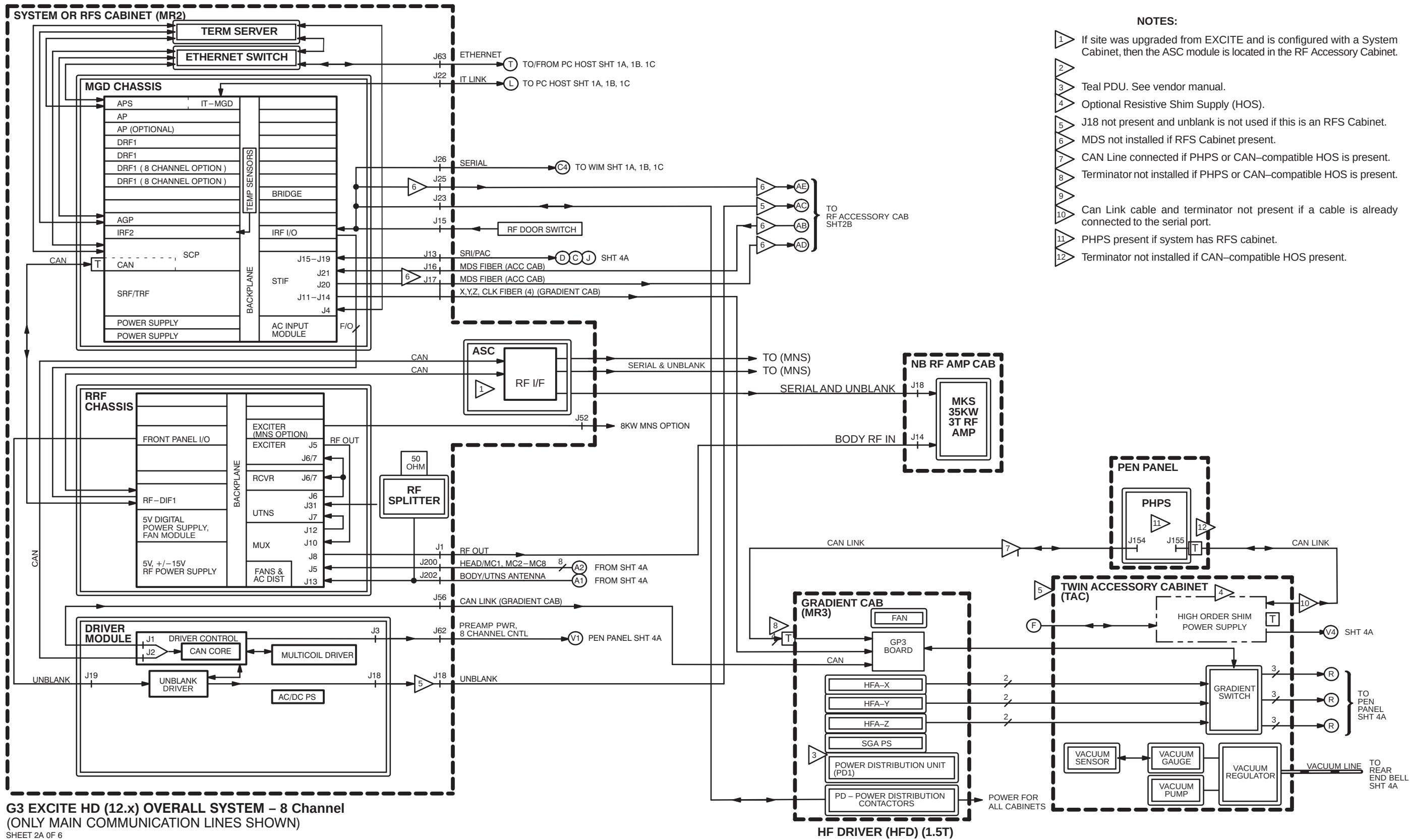



G3 EXCITE HD (12.x) OVERALL SYSTEM - 8 Channel
(ONLY MAIN COMMUNICATION LINES SHOWN)
SHEET 2A OF 6

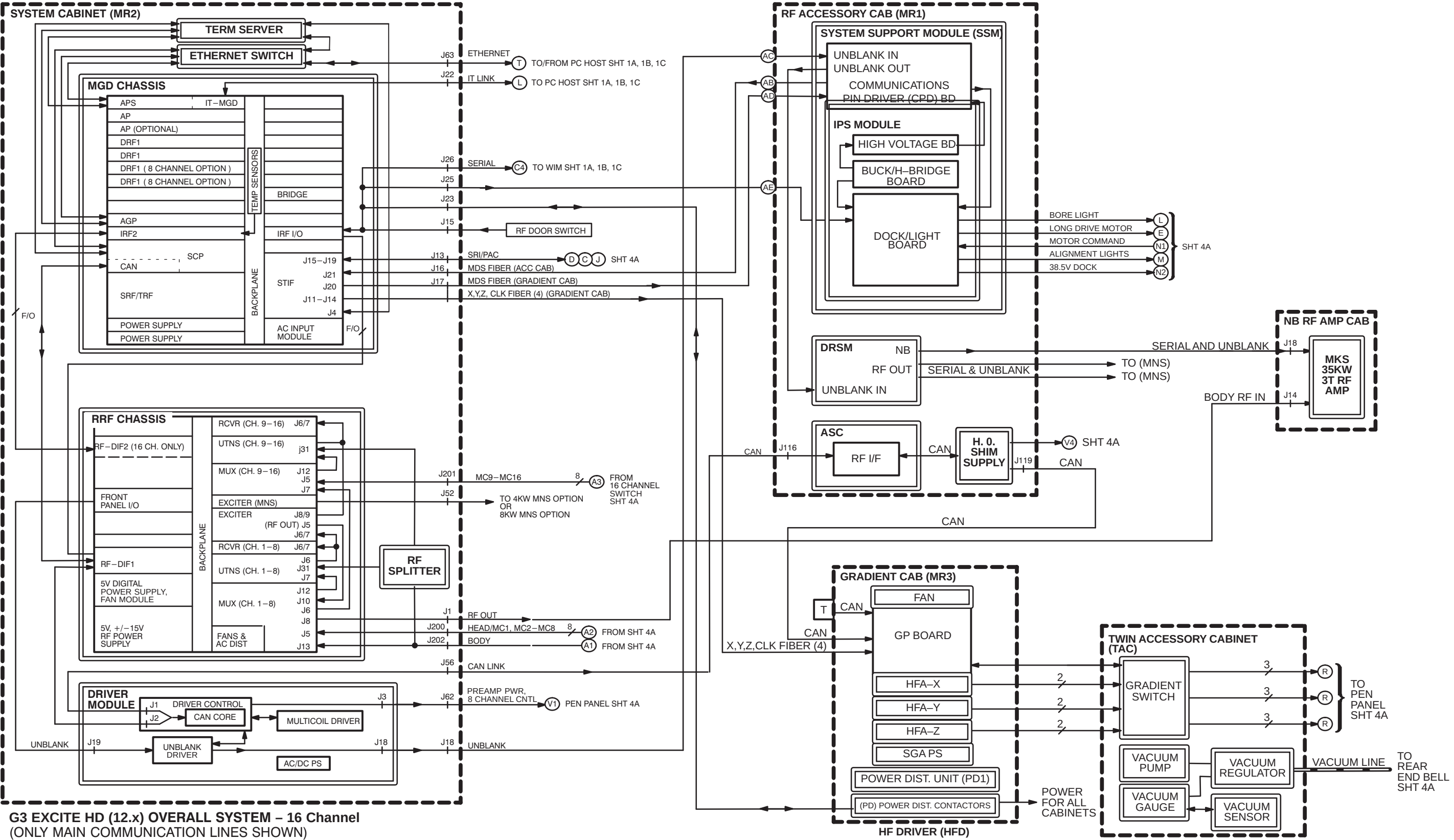

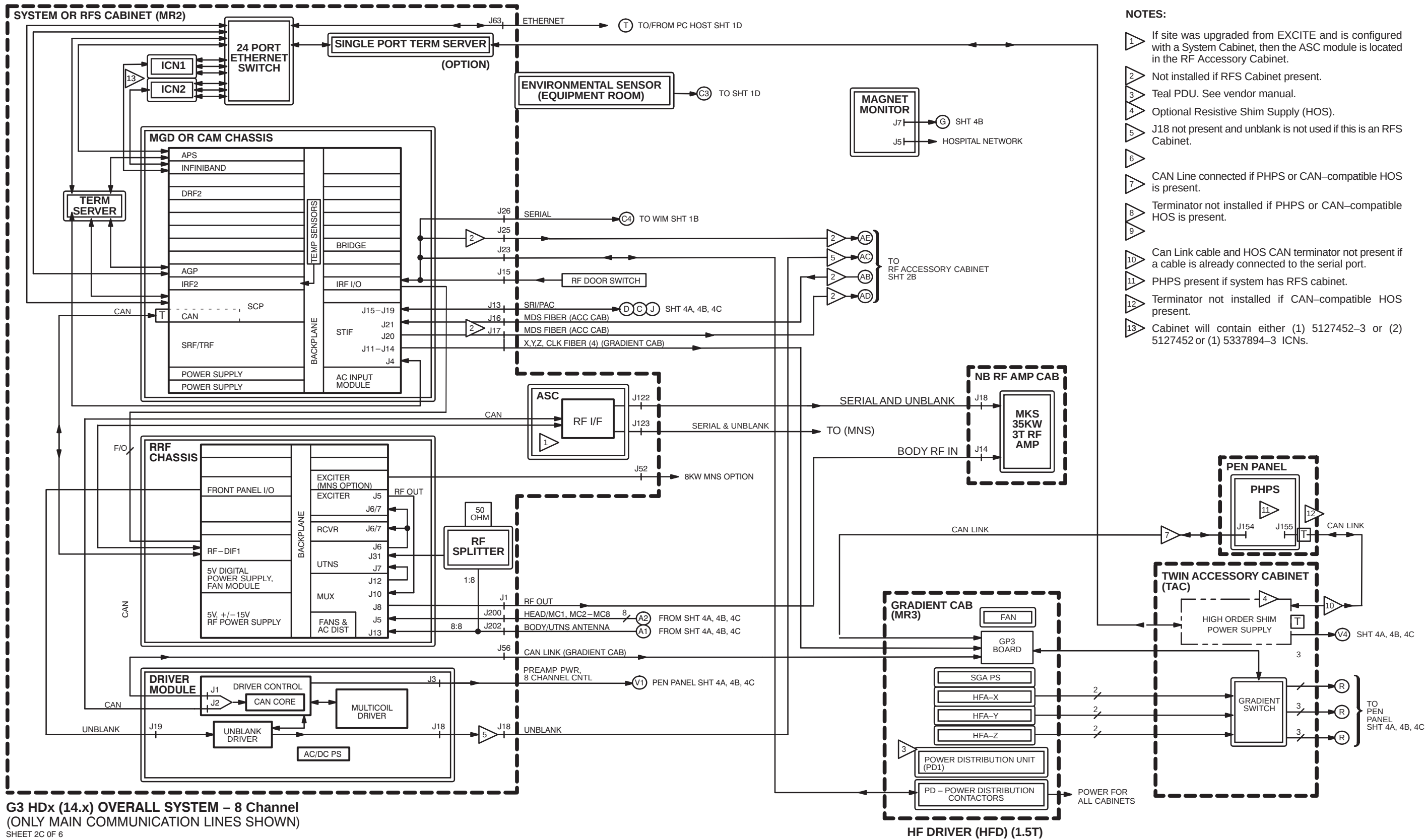



G3 HDx (14.x) OVERALL SYSTEM - 8 Channel
(ONLY MAIN COMMUNICATION LINES SHOWN)
SHEET 2C OF 6

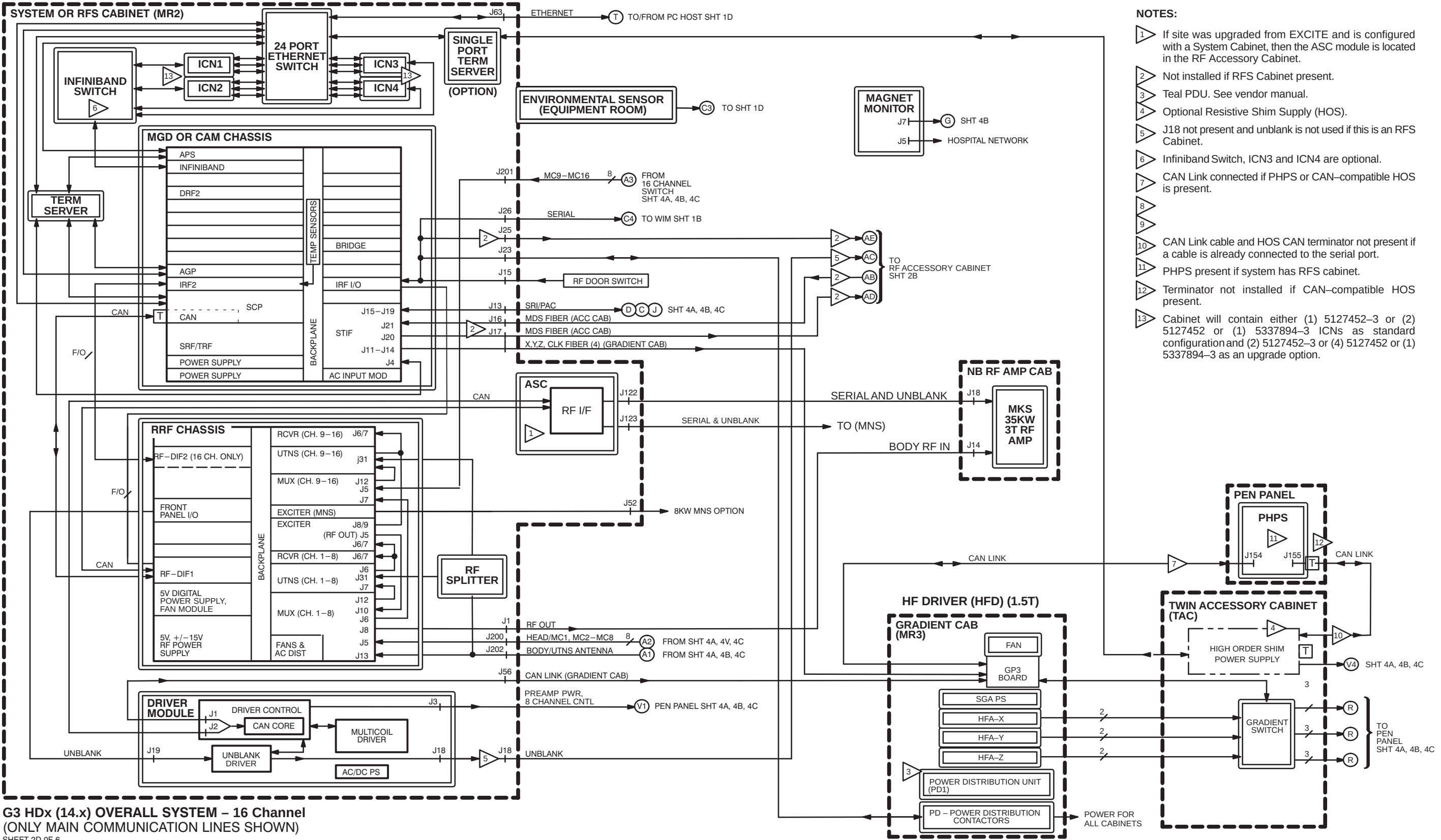

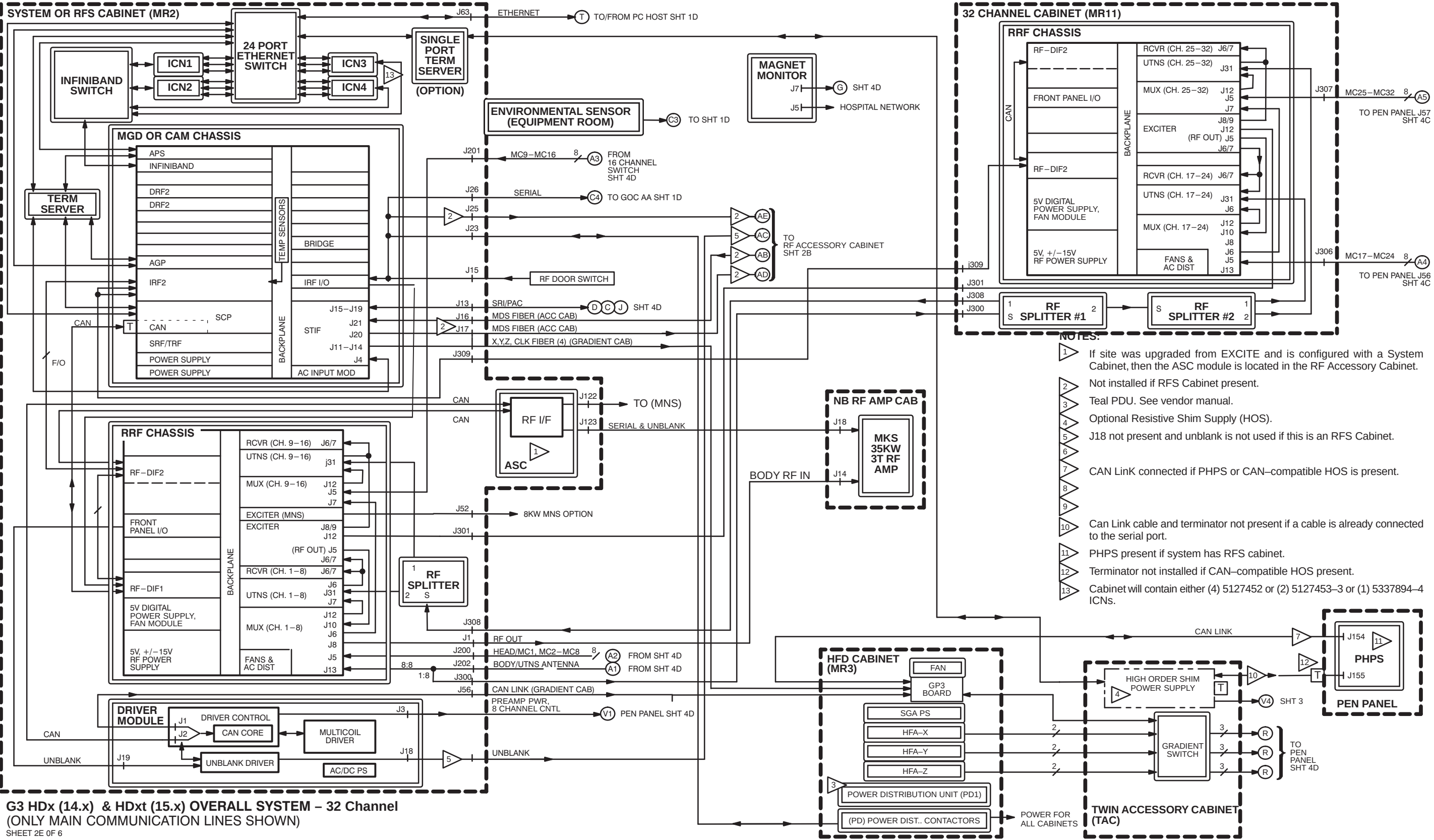


G3 HDx (14.x) & HDxt (15.x) OVERALL SYSTEM – 32 Channel
(ONLY MAIN COMMUNICATION LINES SHOWN)
SHEET 2E OF 6

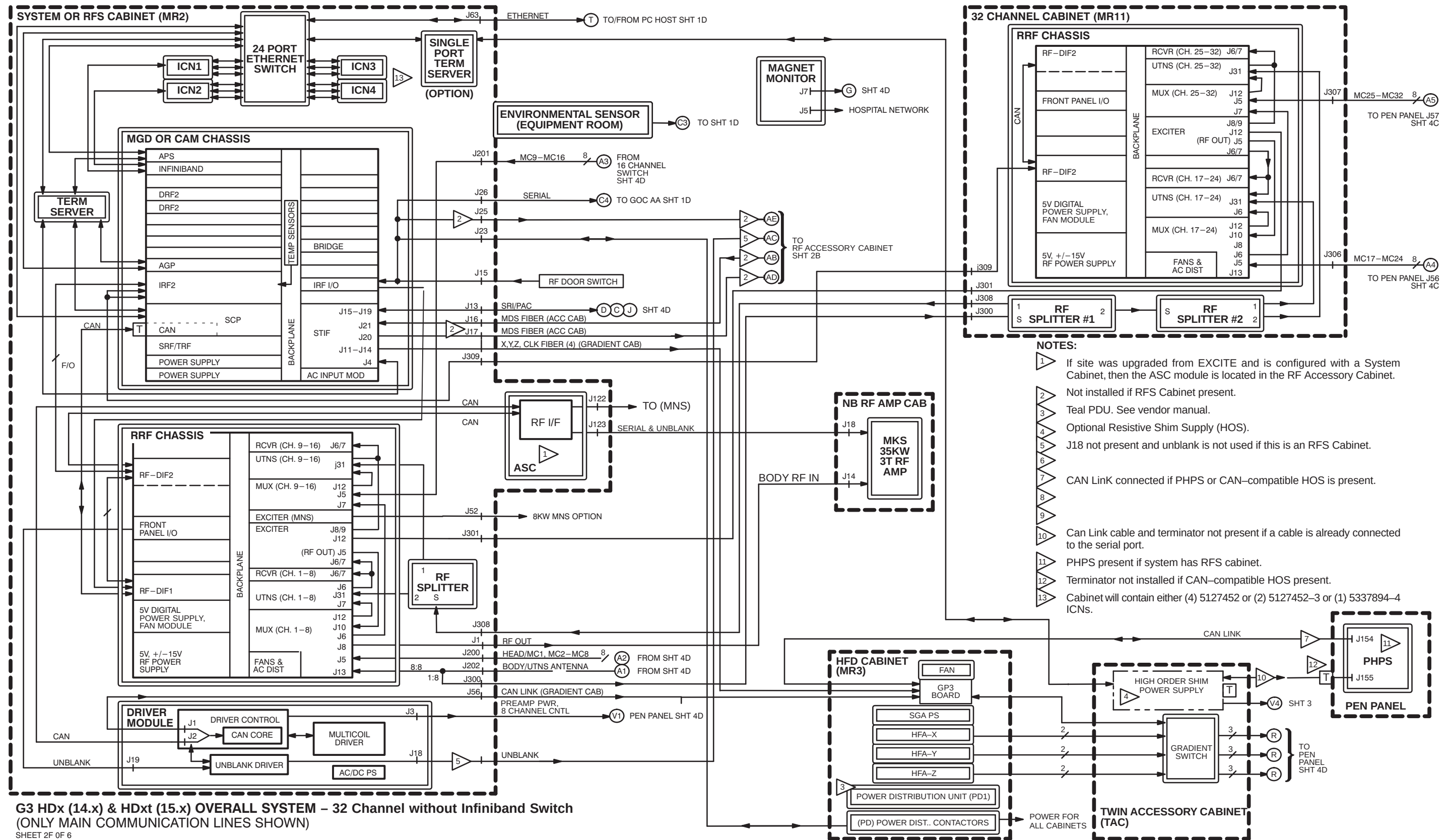

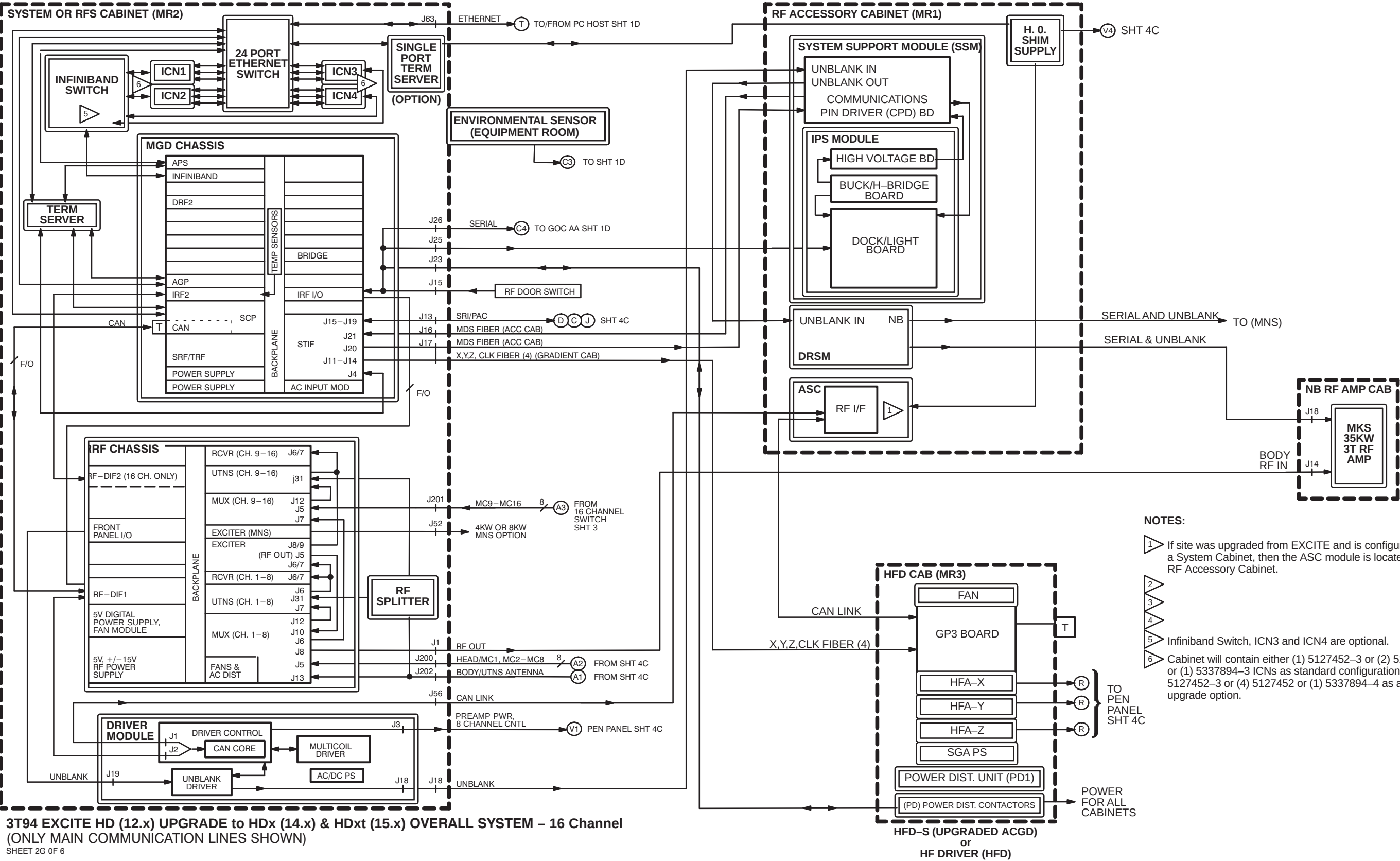



3T94 EXCITE HD (12.x) UPGRADE to HDx (14.x) & HDxt (15.x) OVERALL SYSTEM - 16 Channel
(ONLY MAIN COMMUNICATION LINES SHOWN)


SHEET 2G OF 6

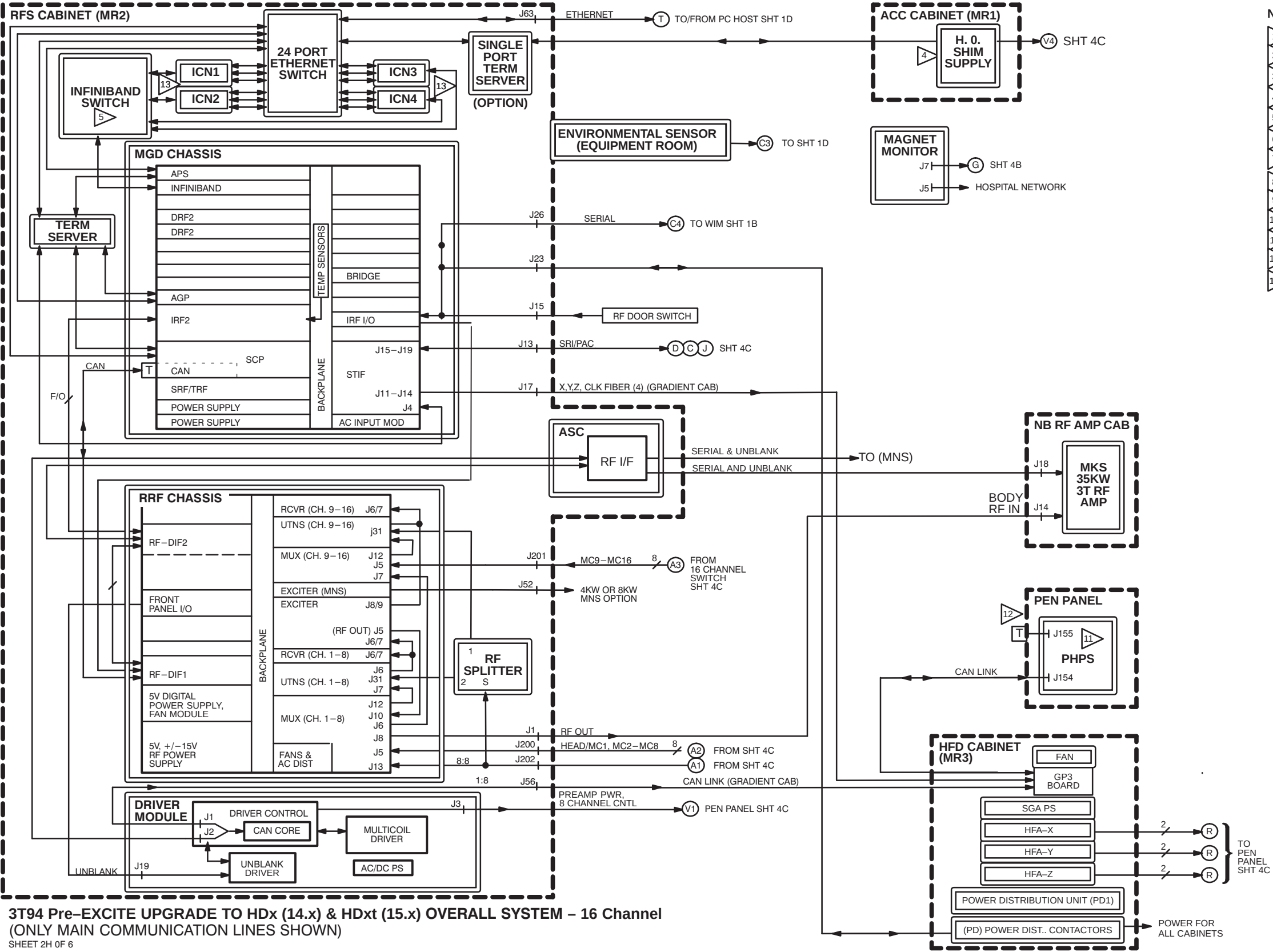

NOTES:

- 1 MDS not installed if RFS Cabinet present.
- 2
- 3
- 4 Optional Resistive Shim Supply (HOS).
- 5 Infiniband Switch, ICN3 and ICN4 are optional.
- 6
- 7
- 8
- 9
- 10
- 11 PPHP present if system has RFS cabinet.
- 12 Terminator not installed if CAN-compatible HOS present.
- 13 Cabinet will contain either (1) 5127452-3 or (2) 5127452 or (1) 5337894-3 as standard configuration and (2) 5127452-3 or (4) 5127452 or (1) 5337894-3 as an upgrade option.

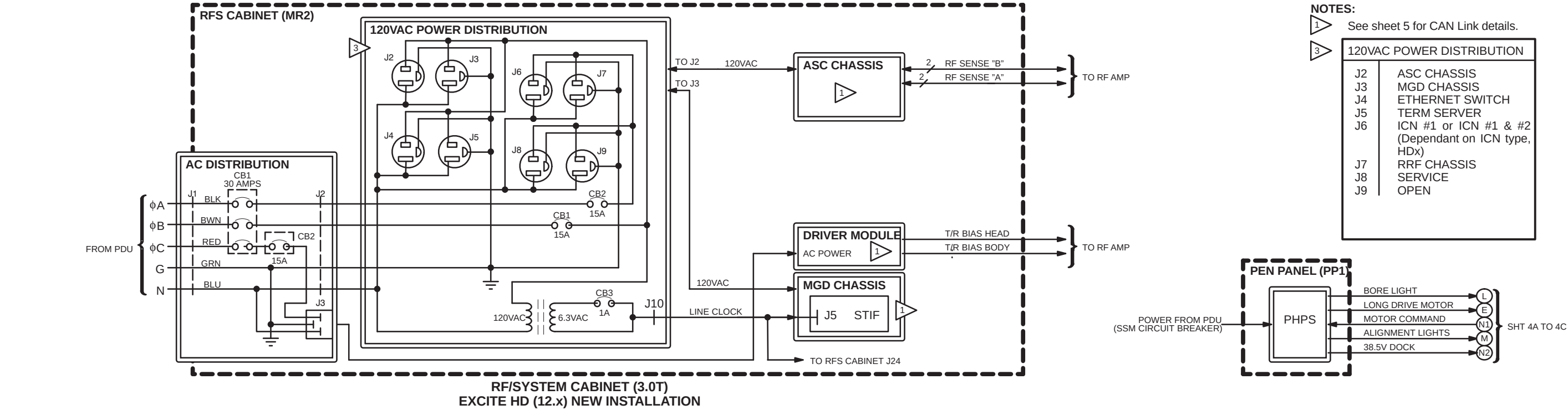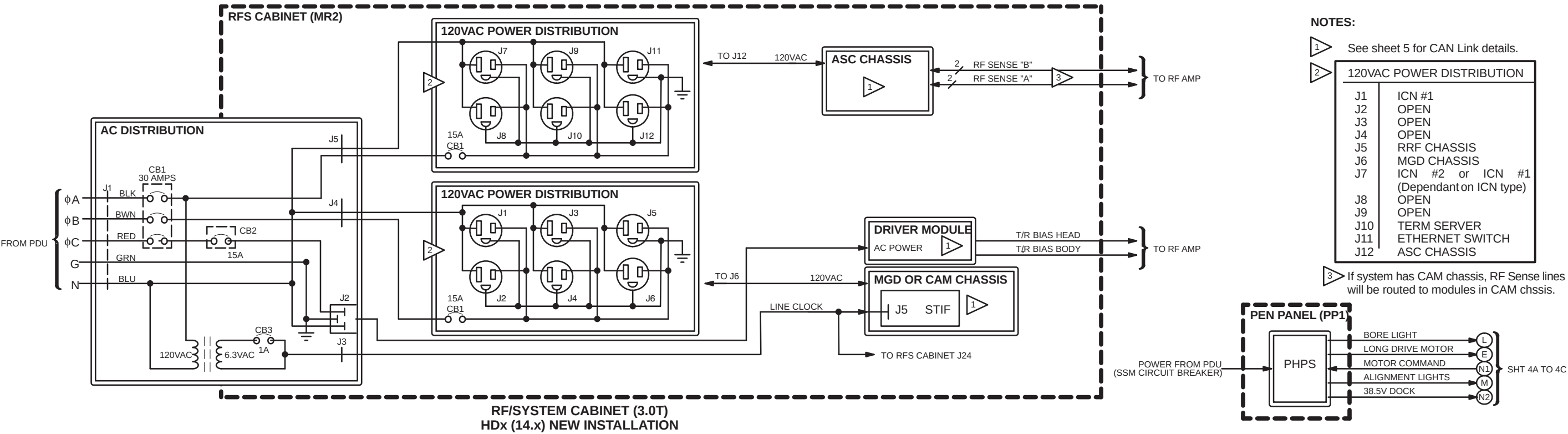

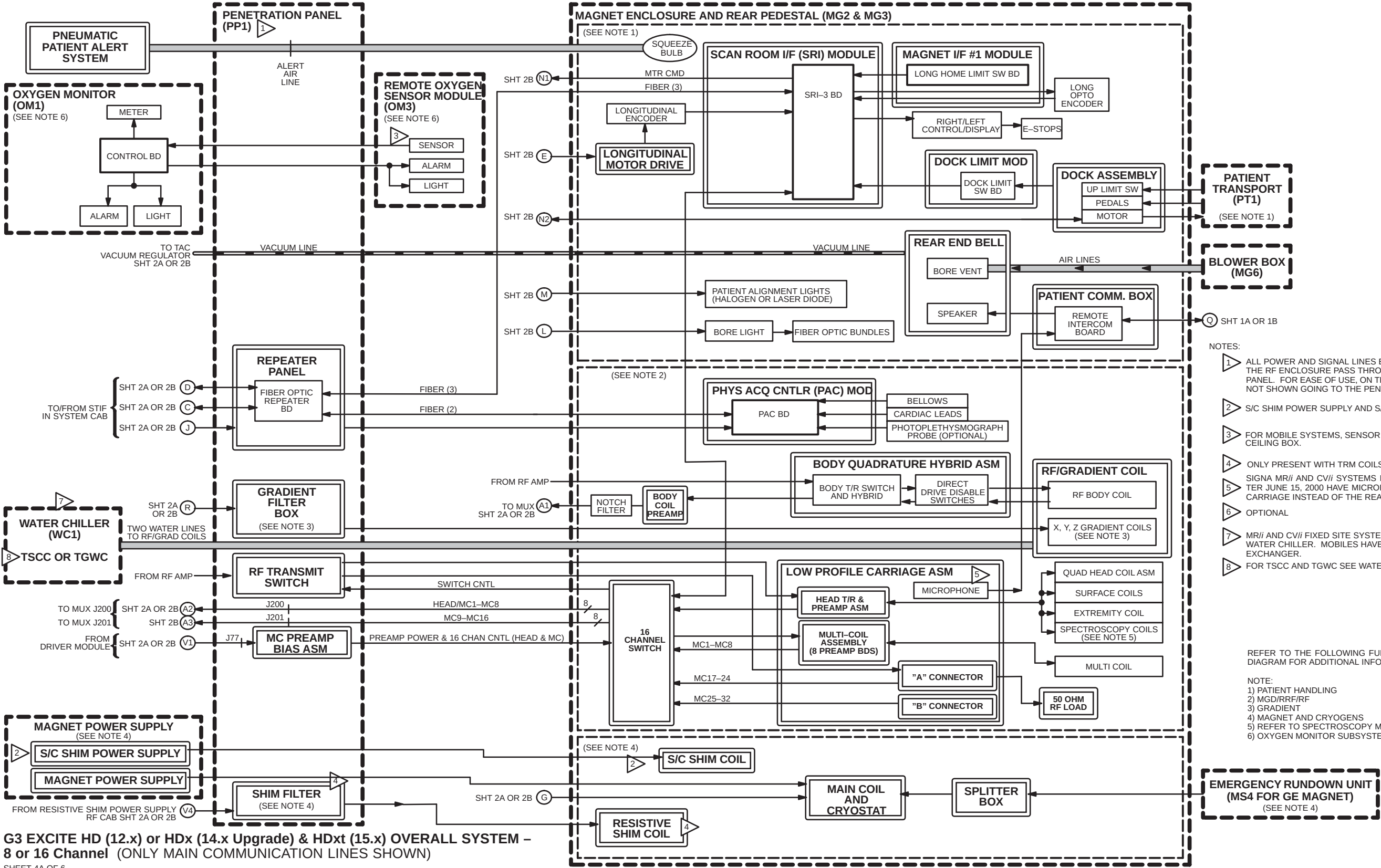


G3 EXCITE HD (12.x) or HDx (14.x Upgrade) & HDxt (15.x) OVERALL SYSTEM –
8 or 16 Channel (ONLY MAIN COMMUNICATION LINES SHOWN)

SHEET 4A OF 6


- NOTES:
- 1 ALL POWER AND SIGNAL LINES ENTERING AND LEAVING THE RF ENCLOSURE PASS THROUGH THE PENETRATION PANEL. FOR EASE OF USE, ON THIS DIAGRAM, THEY ARE NOT SHOWN GOING TO THE PENETRATION PANEL.
 - 2 S/C SHIM POWER SUPPLY AND S/C SHIM COILS.
 - 3 FOR MOBILE SYSTEMS, SENSOR IS LOCATED IN CEILING BOX.
 - 4 ONLY PRESENT WITH TRM COILS.
 - 5 SIGNA MR*ii* AND CV*ii* SYSTEMS MANUFACTURED AFTER JUNE 15, 2000 HAVE MICROPHONE IN THE HEAD CARRIAGE INSTEAD OF THE REAR END BELL.
 - 6 OPTIONAL
 - 7 MR*ii* AND CV*ii* FIXED SITE SYSTEMS HAVE NESLAB WATER CHILLER. MOBILES HAVE E & W HEAT EXCHANGER.
 - 8 FOR TSCC AND TGWC SEE WATER CHILLER.

REFER TO THE FOLLOWING FUNCTIONAL BLOCK DIAGRAM FOR ADDITIONAL INFORMATION:

- NOTE:
- 1) PATIENT HANDLING
 - 2) MGD/RRF/RF
 - 3) GRADIENT
 - 4) MAGNET AND CRYOGENS
 - 5) REFER TO SPECTROSCOPY MANUALS
 - 6) OXYGEN MONITOR SUBSYSTEM (DIR. 15336).

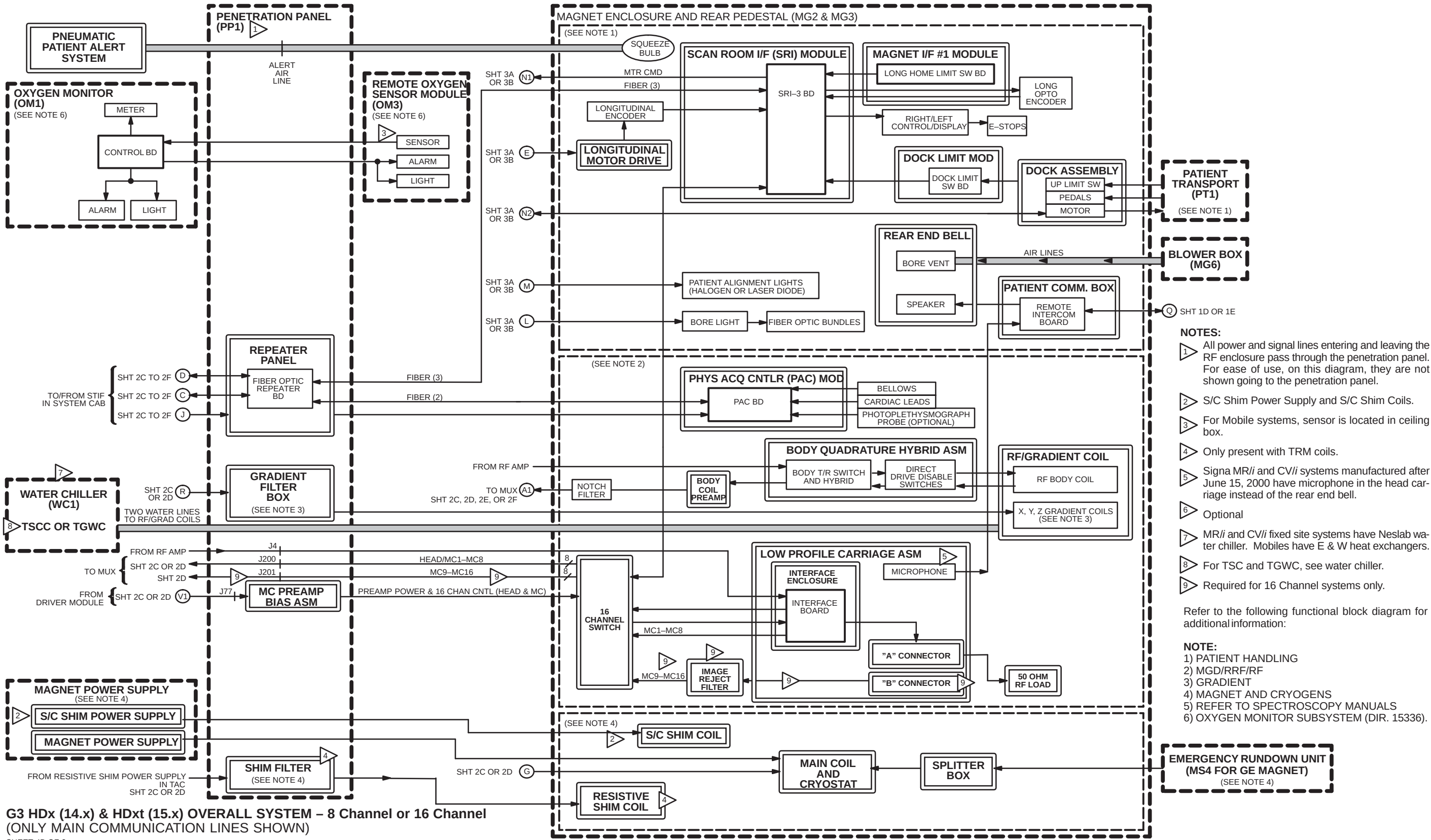


G3 HDx (14.x) & HDxt (15.x) OVERALL SYSTEM – 8 Channel or 16 Channel
(ONLY MAIN COMMUNICATION LINES SHOWN)
SHEET 4B OF 6

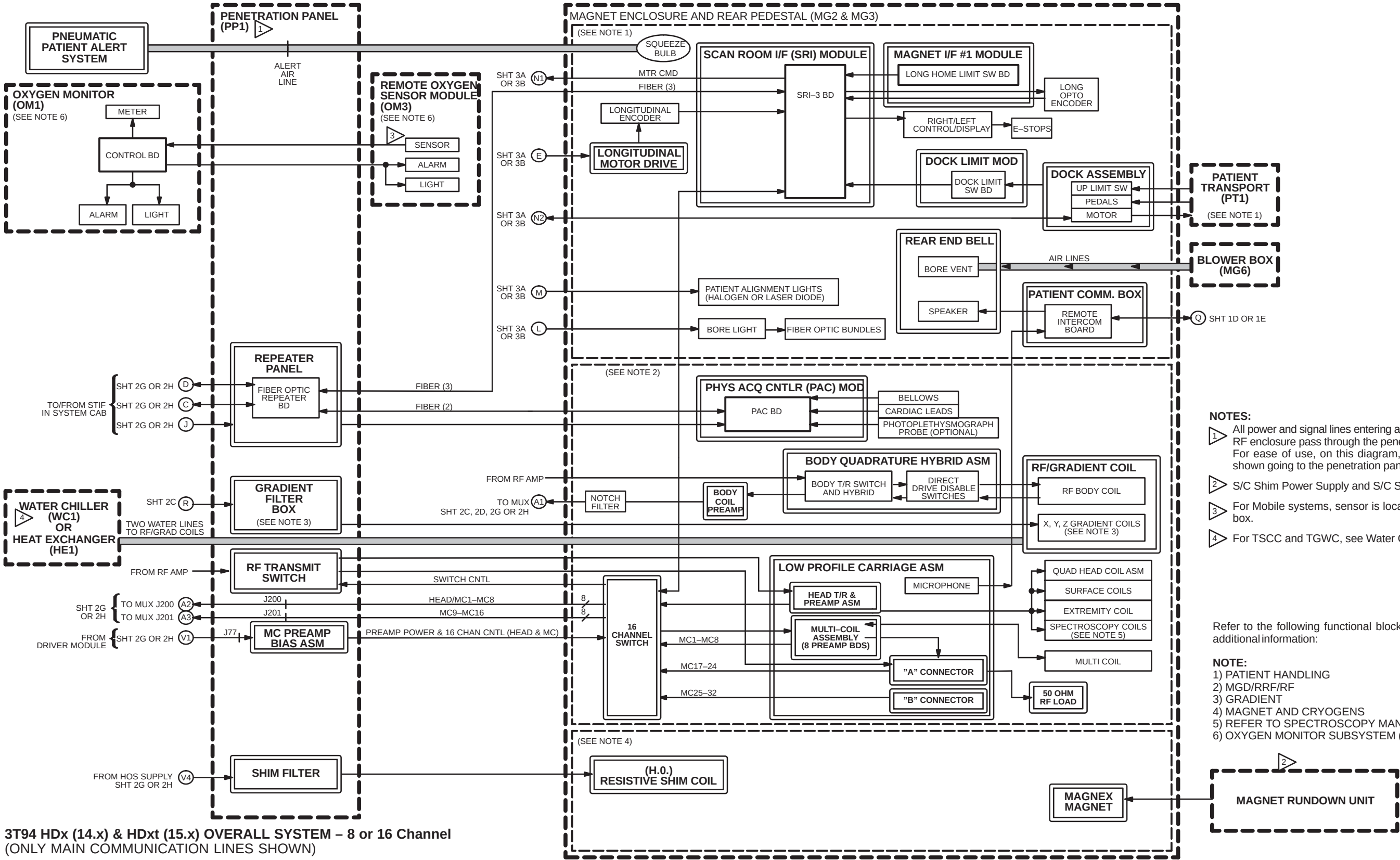


3T94 HDx (14.x) & HDxt (15.x) OVERALL SYSTEM – 8 or 16 Channel
(ONLY MAIN COMMUNICATION LINES SHOWN)
SHEET 4C OF 6

- NOTES:**
- 1 All power and signal lines entering and leaving the RF enclosure pass through the penetration panel. For ease of use, on this diagram, they are not shown going to the penetration panel.
 - 2 S/C Shim Power Supply and S/C Shim Coils.
 - 3 For Mobile systems, sensor is located in ceiling box.
 - 4 For TSCC and TGWC, see Water Chiller.


- Refer to the following functional block diagram for additional information:
- NOTE:**
- 1) PATIENT HANDLING
 - 2) MGD/RRF/RF
 - 3) GRADIENT
 - 4) MAGNET AND CRYOGENS
 - 5) REFER TO SPECTROSCOPY MANUALS
 - 6) OXYGEN MONITOR SUBSYSTEM (DIR. 15336).

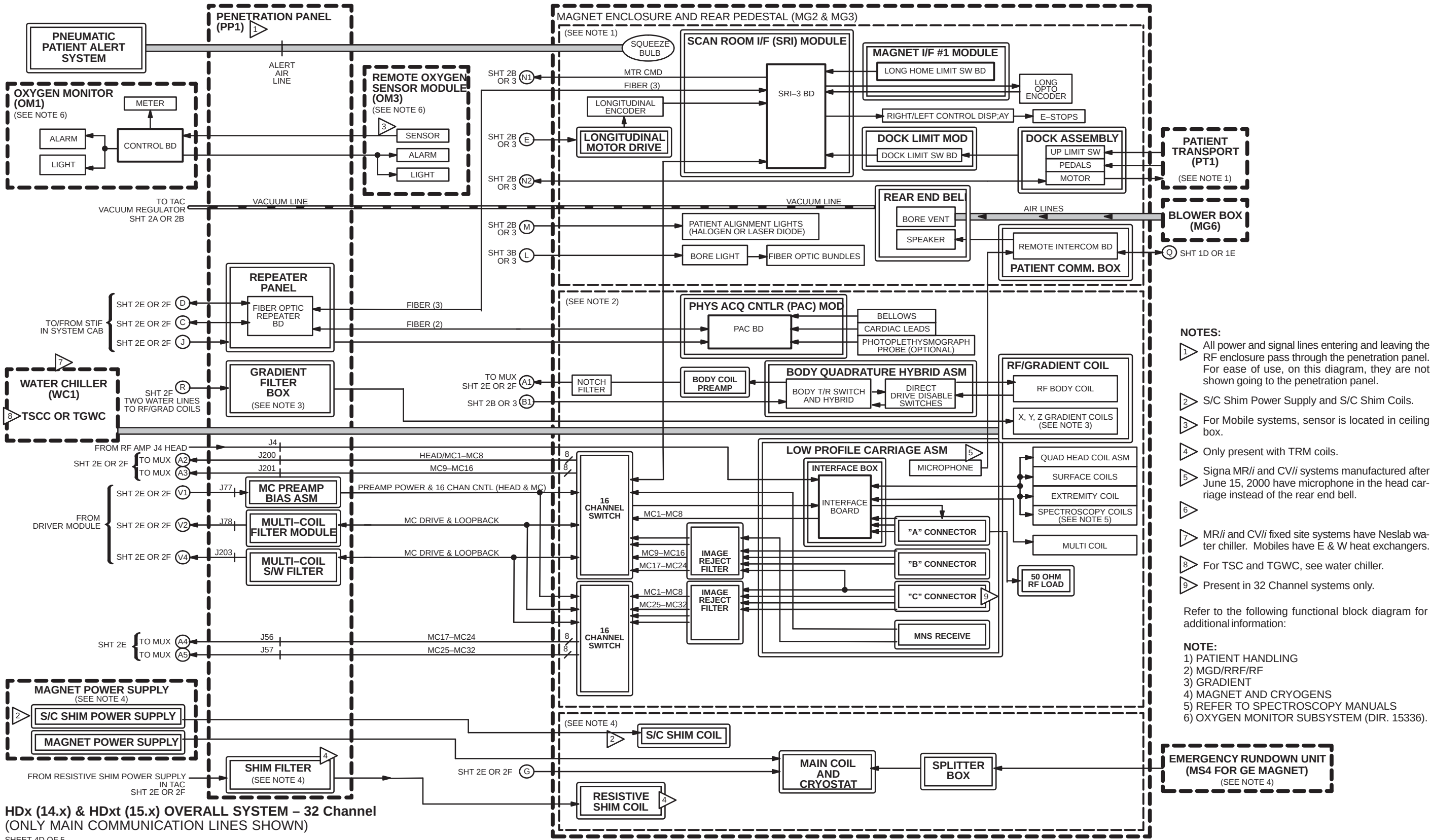


HDx (14.x) & HDxt (15.x) OVERALL SYSTEM – 32 Channel
(ONLY MAIN COMMUNICATION LINES SHOWN)
SHEET 4D OF 5

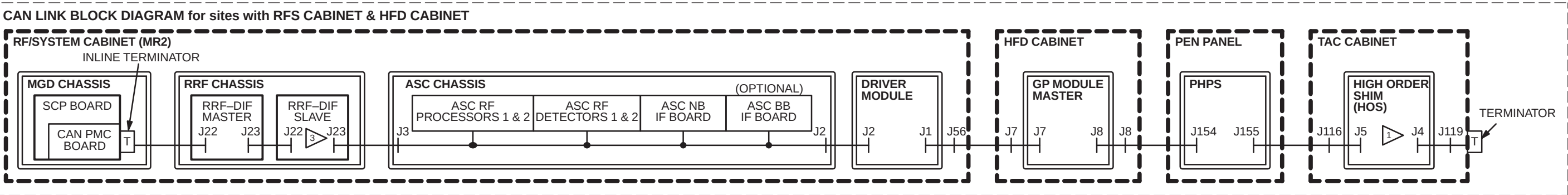

OR

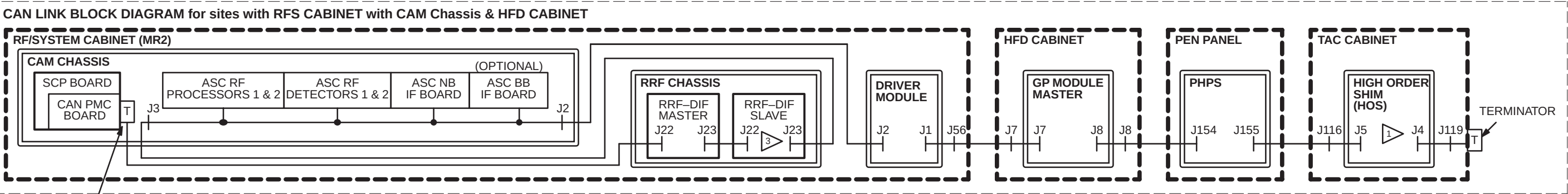


OR

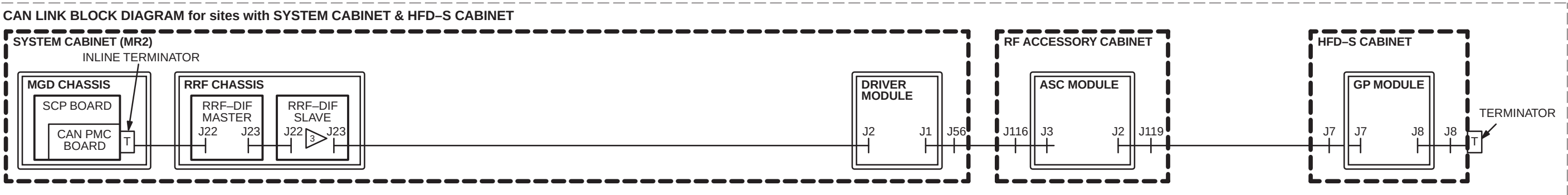

- NOTES:
- 1

HOS is optional. Terminator must be connected to J155 on PHPS if HOS not present.
- 2
- 3

Only present in systems equipped with 16 or more Channels. If system is 8 Channels or less, then RF-DIF Master J23 connects to next unit in link as shown.

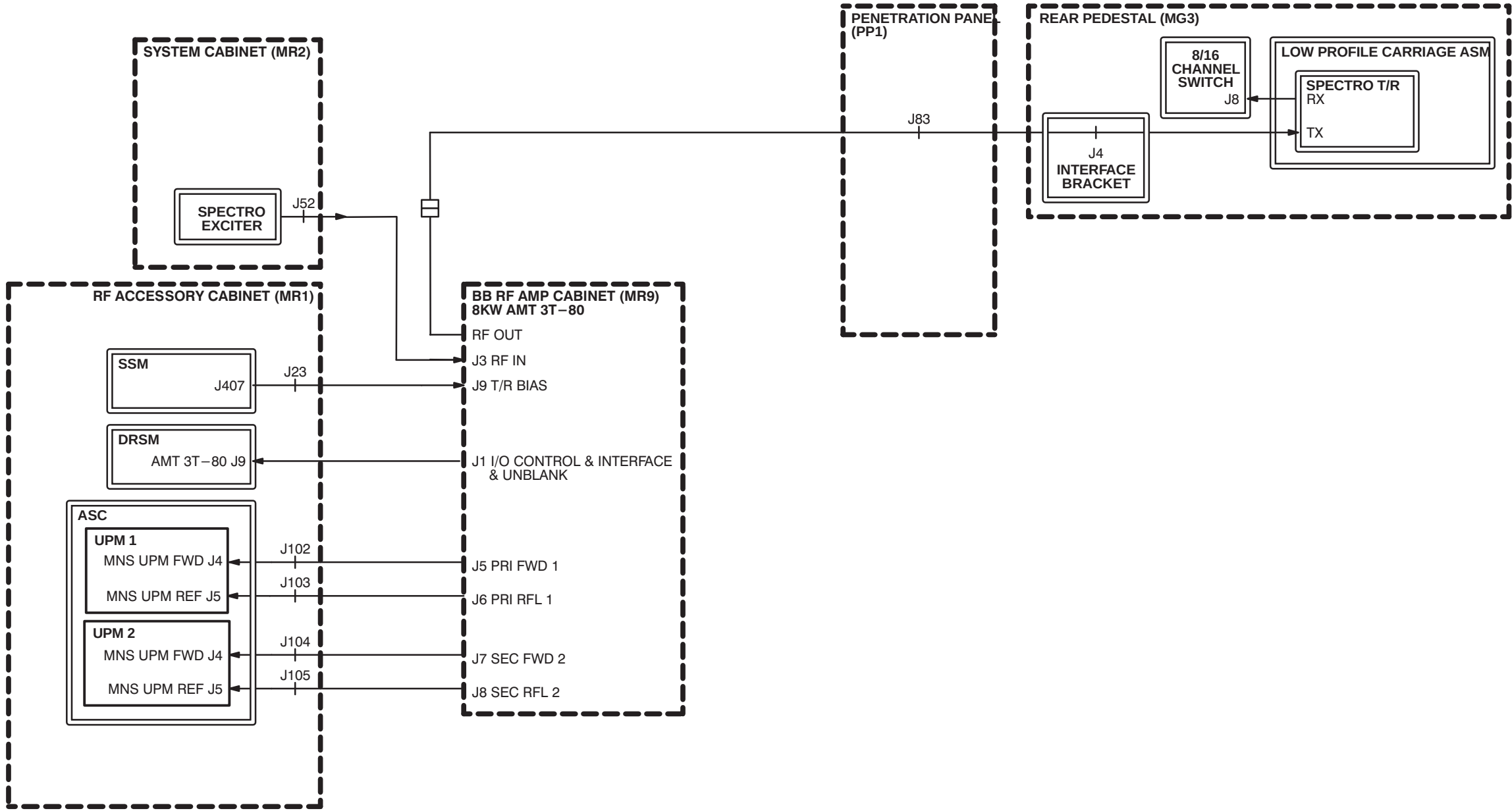


G3 EXCITE HD (12.x) and HDx (14.x) & HDxt (15.x) OVERALL SYSTEM – 8KW MNS OPTION
(ONLY MAIN COMMUNICATION LINES SHOWN)
SHEET 5 OF 6


OPERATOR WORKSPACE

TABLE OF CONTENTS

<u>SECTION</u>	<u>PAGE</u>
1 – SYSTEM	
G3 & 3T94 HD OVERALL SYSTEM – OPERATOR WORKSPACE CABINET w/8000 & GOC AA (Forward Production or Pre–EXCITE to HD Upgrades)	1–1
G3 & 3T94 HD OVERALL SYSTEM – OPERATOR WORKSPACE CABINET w/8000 & WIM (EXCITE to HD Upgrades)	1–2
G3 & 3T94 HD OVERALL SYSTEM – OPERATOR WORKSPACE CABINET w/8200 & GOC AA (Forward Production or Pre–EXCITE to HD Upgrades)	1–3
G3 & 3T94 HDx & HDxt OVERALL SYSTEM – OPERATOR WORKSPACE CABINETw/9300 & GOC AA (Forward Production, EXCITE, or Pre–EXCITE to HDx Upgrades)	1–4
G3 & 3T94 HD OVERALL SYSTEM – OPERATOR WORKSPACE CABINET/TABLE w/WIM (Forward Production & Upgrades)	1–5
G3 & 3T94 HD OVERALL SYSTEM – OPERATOR WORKSPACE CABINET/TABLE w/GOC AA (Forward Production & Upgrades)	1–6
G3 & 3T94 HDx & HDxt OVERALL SYSTEM – OPERATOR WORKSPACE CABINET/TABLE w/GOC AA (Forward Production & Upgrades)	1–7

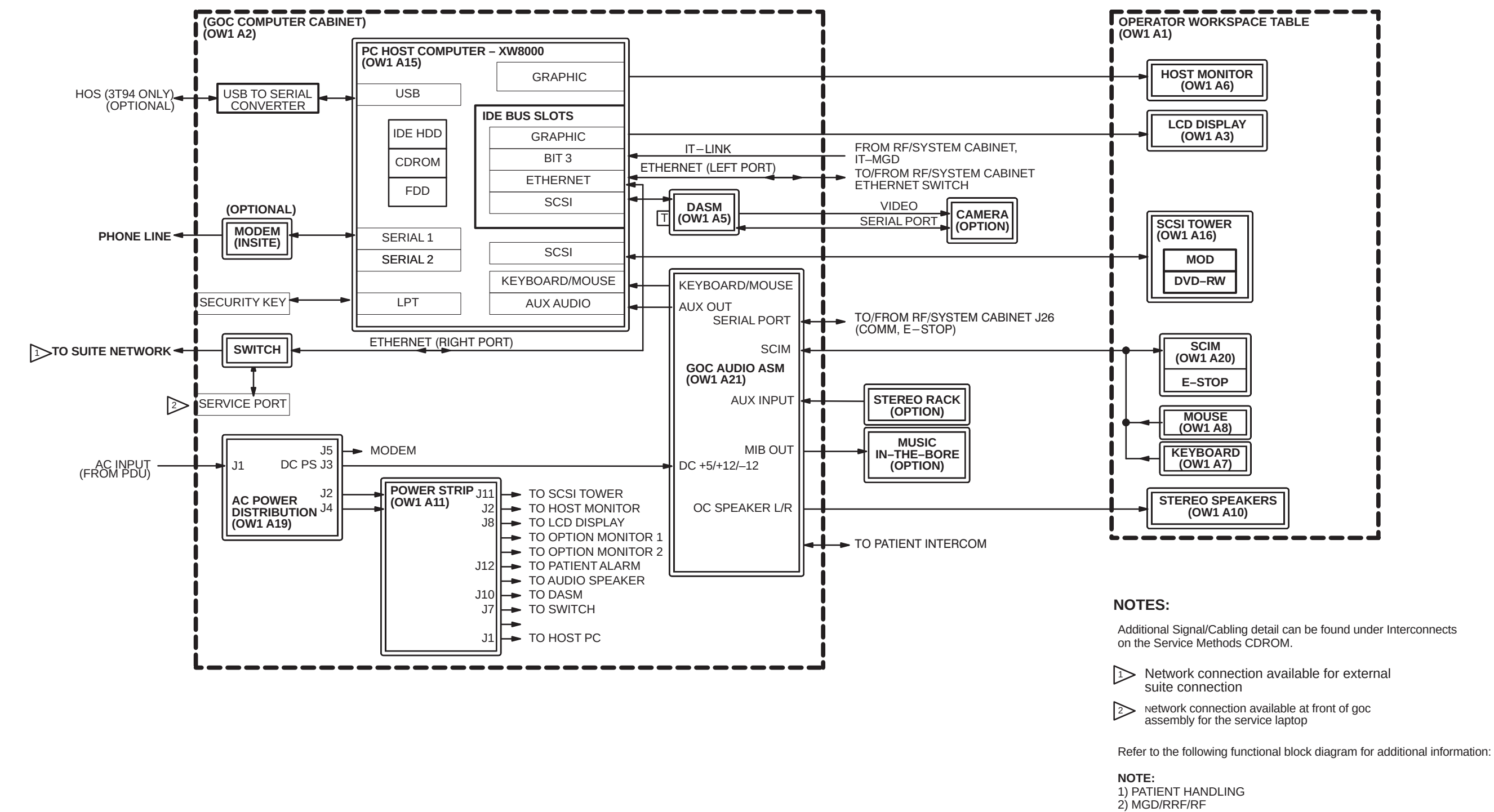

G3 & 3T94 EXCITE HD (12.x) OVERALL SYSTEM
Forward Production or Pre-EXCITE Upgrades
(OPERATOR WORKSPACE CABINET w/8000 & GOC AA and TABLE)
SHEET 1A OF 5

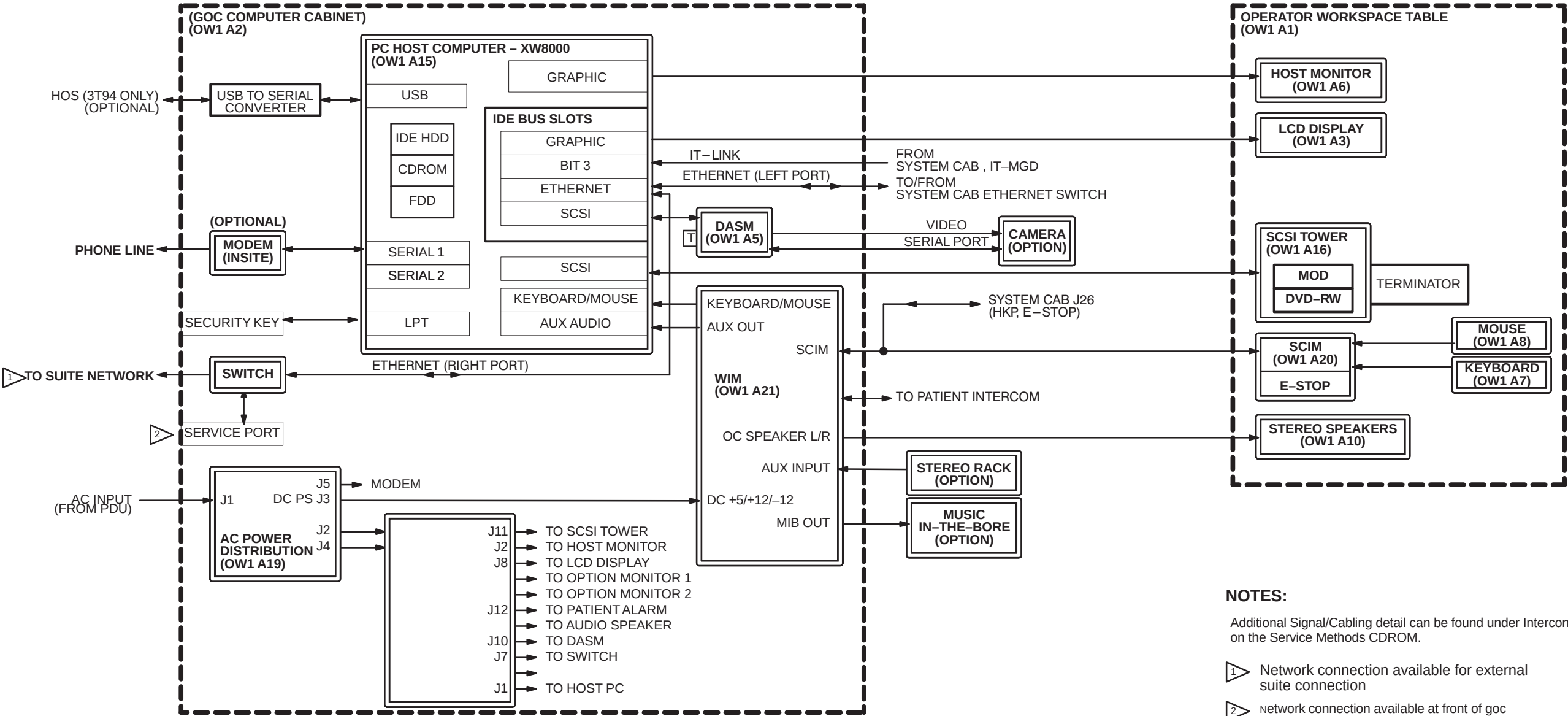

NOTES:

Additional Signal/Cabling detail can be found under Interconnects on the Service Methods CDROM.

1 Network connection available for external suite connection


2 network connection available at front of goc assembly for the service laptop

Refer to the following functional block diagram for additional information:


NOTE:

1) PATIENT HANDLING

2) MGD/RRF/RF

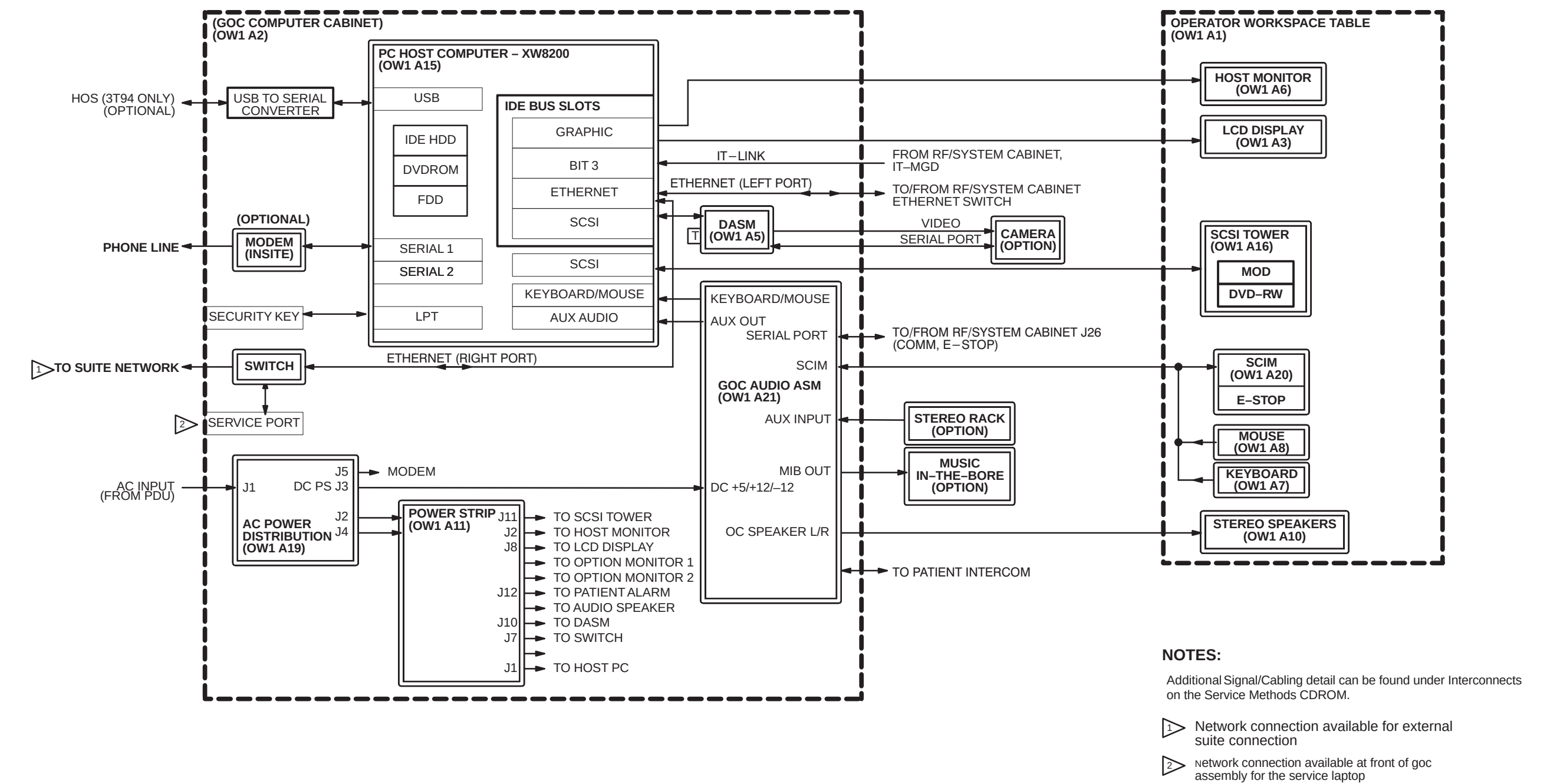

G3 & 3T94 EXCITE HD (12.x) OVERALL SYSTEM
Forward Production or Pre-EXCITE Upgrades
(OPERATOR WORKSPACE CABINET w/8200 & GOC AA and TABLE)
SHEET 1C OF 5

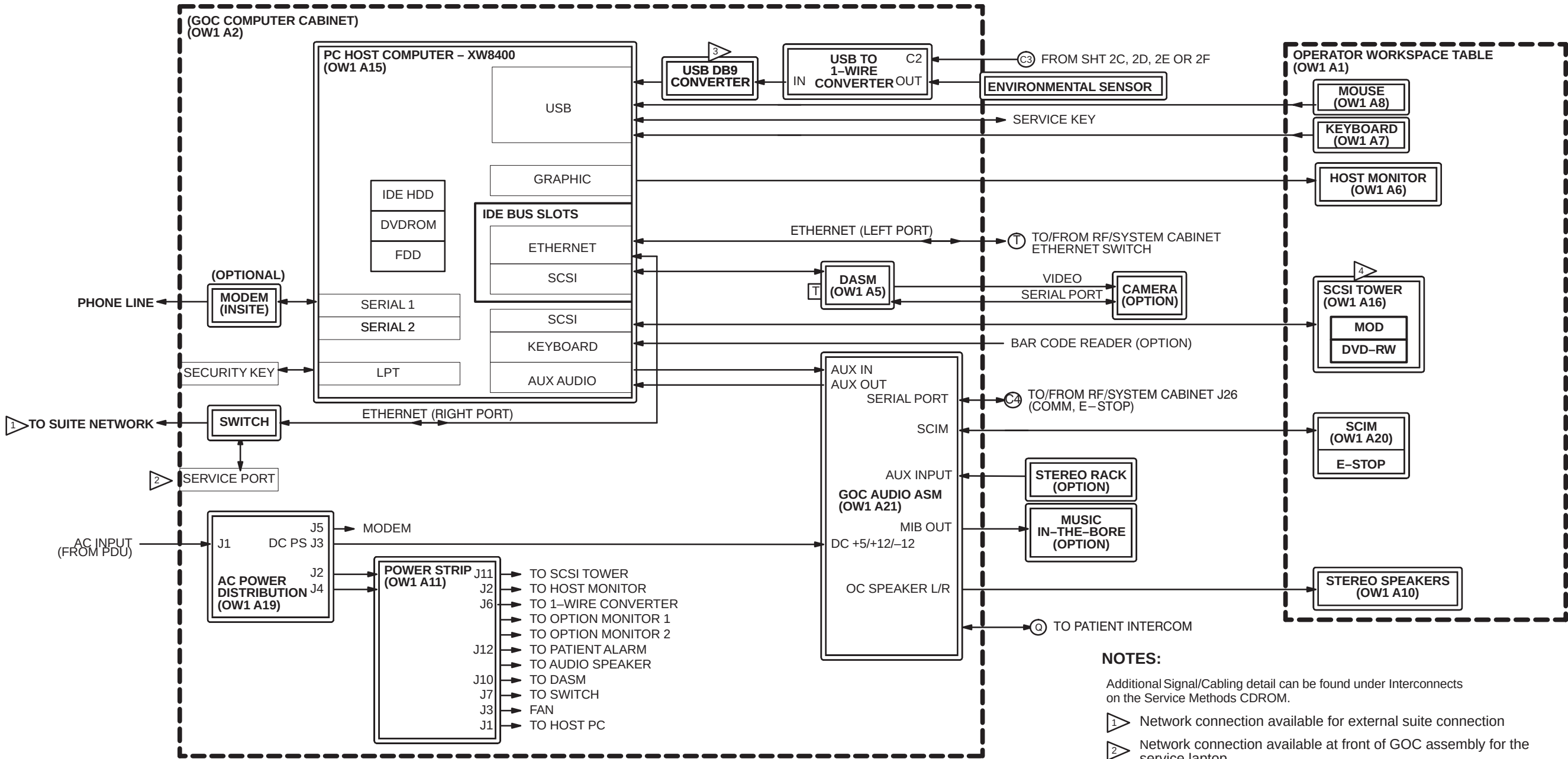


Refer to the following functional block diagram for additional information:

NOTE:
1) PATIENT HANDLING
2) MGD/RRF/RF


NOTES:

Additional Signal/Cabling detail can be found under Interconnects on the Service Methods CDROM.

- 1 Network connection available for external suite connection
- 2 Network connection available at front of GOC assembly for the service laptop
- 3 Beginning August 2008, design changed to 1-wire to USB converter. Environmental Sensors Block Diagram utilizing 1-wire 5147136-2
- 4 8400 Host Computer is also available without SCSI Tower.

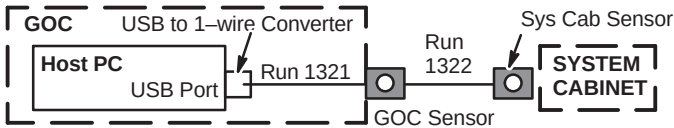

HDx (14.x) & HDxt (15.x) OVERALL SYSTEM
Forward Production, EXCITE, or Pre-EXCITE to HDx Upgrades
(OPERATOR WORKSPACE CABINET w/8400 & TABLE)

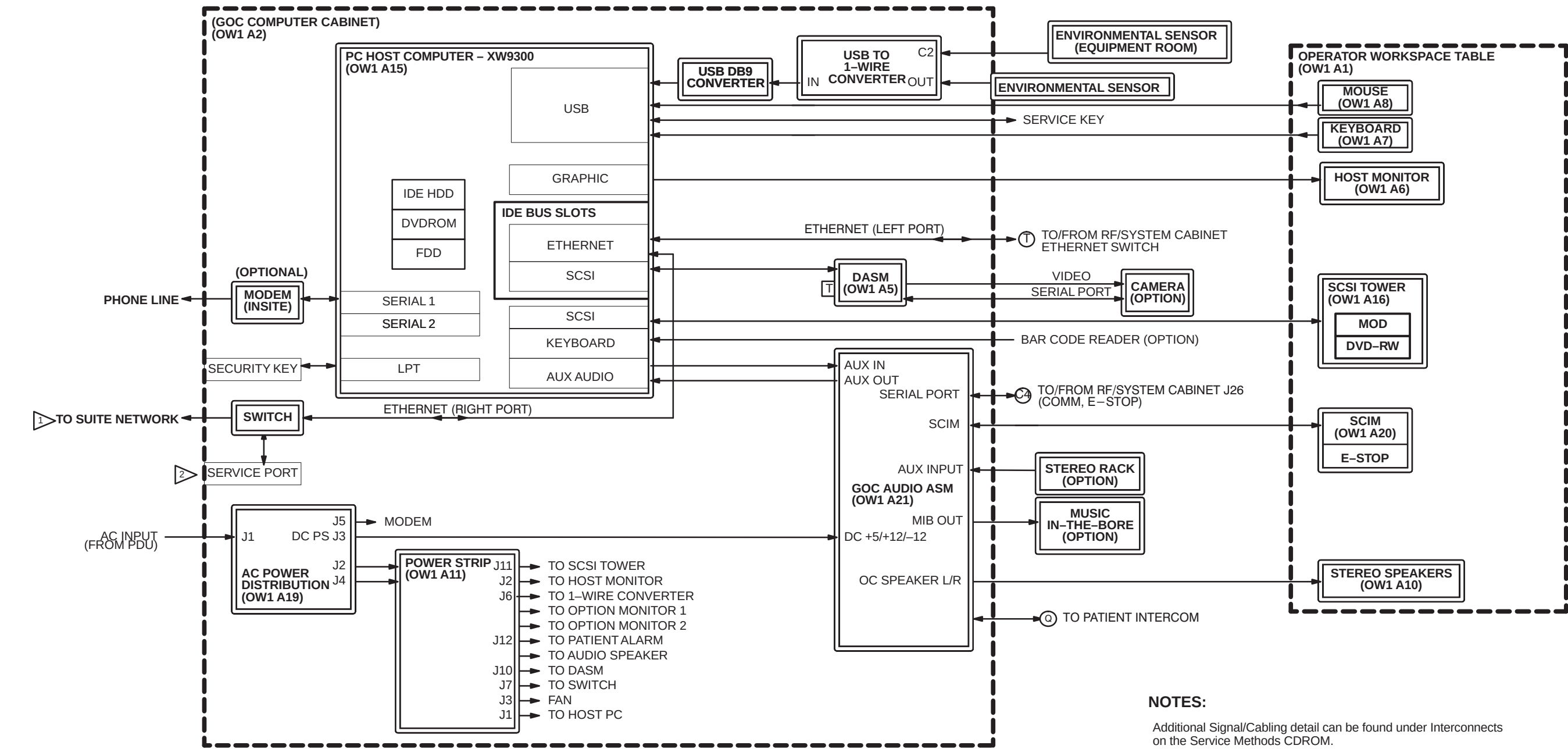

NOTES:

Additional Signal/Cabling detail can be found under Interconnects on the Service Methods CDROM.

1 Network connection available for external suite connection


2 Network connection available at front of GOC assembly for the service laptop

Refer to the following functional block diagram for additional information:

NOTE:

1) PATIENT HANDLING

2) MGD/RRF/RF

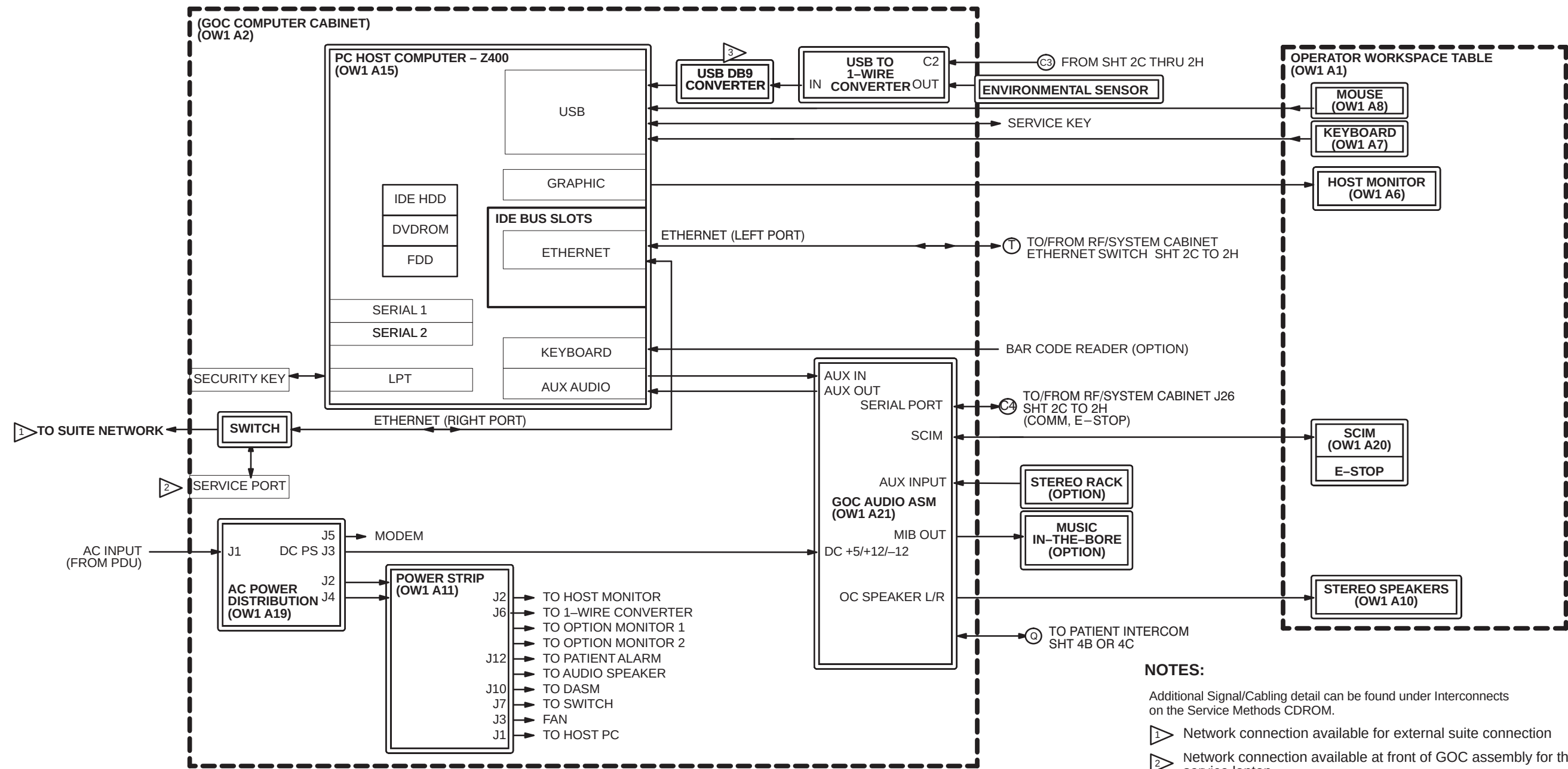

NOTES:

Additional Signal/Cabling detail can be found under Interconnects on the Service Methods CDROM.


- 1 Network connection available for external suite connection
- 2 Network connection available at front of GOC assembly for the service laptop
- 3 Beginning August 2008, design changed to 1-wire to USB converter. Environmental Sensors Block Diagram utilizing 1-wire 5147136-2

The small schematic shows a **Host PC** with a **USB Port** connected to a **USB to 1-wire Converter**. This converter is connected to **Run 1321**, which leads to a **GOC Sensor**. **Run 1322** connects the **GOC Sensor** to a **Sys Cab Sensor** located in the **SYSTEM CABINET**.


G3 & 3T94 HDx (14.x) & HDxt (15.x) OVERALL SYSTEM
Forward Production, EXCITE, or Pre-EXCITE to HDx Upgrades
(OPERATOR WORKSPACE CABINET w/Z400 & TABLE)
SHEET 1F OF 6

Refer to the following functional block diagram for additional information:

NOTE:
1) PATIENT HANDLING
2) MGD/RRF/RF

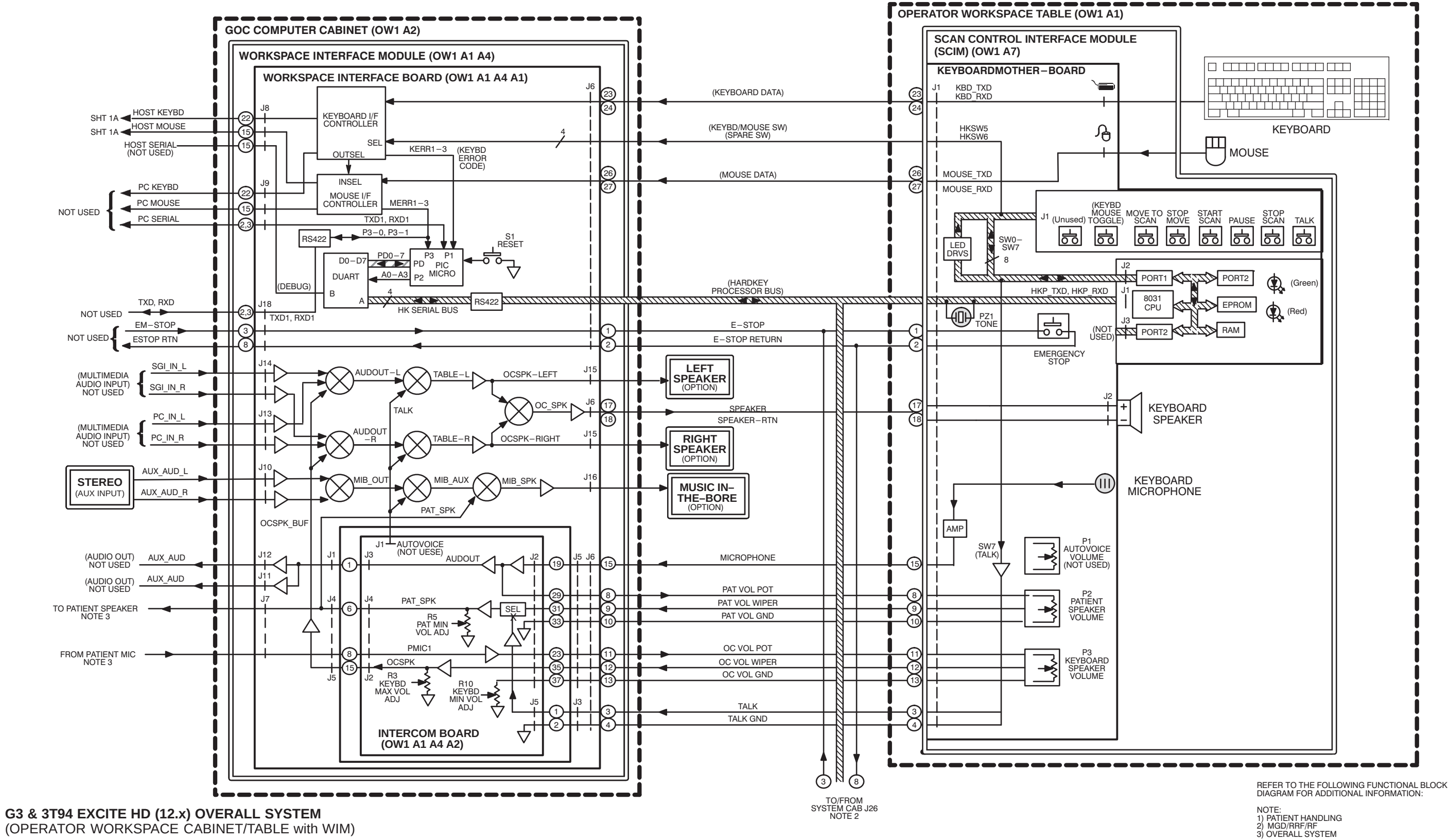

G3 & 3T94 EXCITE HD (12.x) OVERALL SYSTEM
(OPERATOR WORKSPACE CABINET/TABLE with WIM)

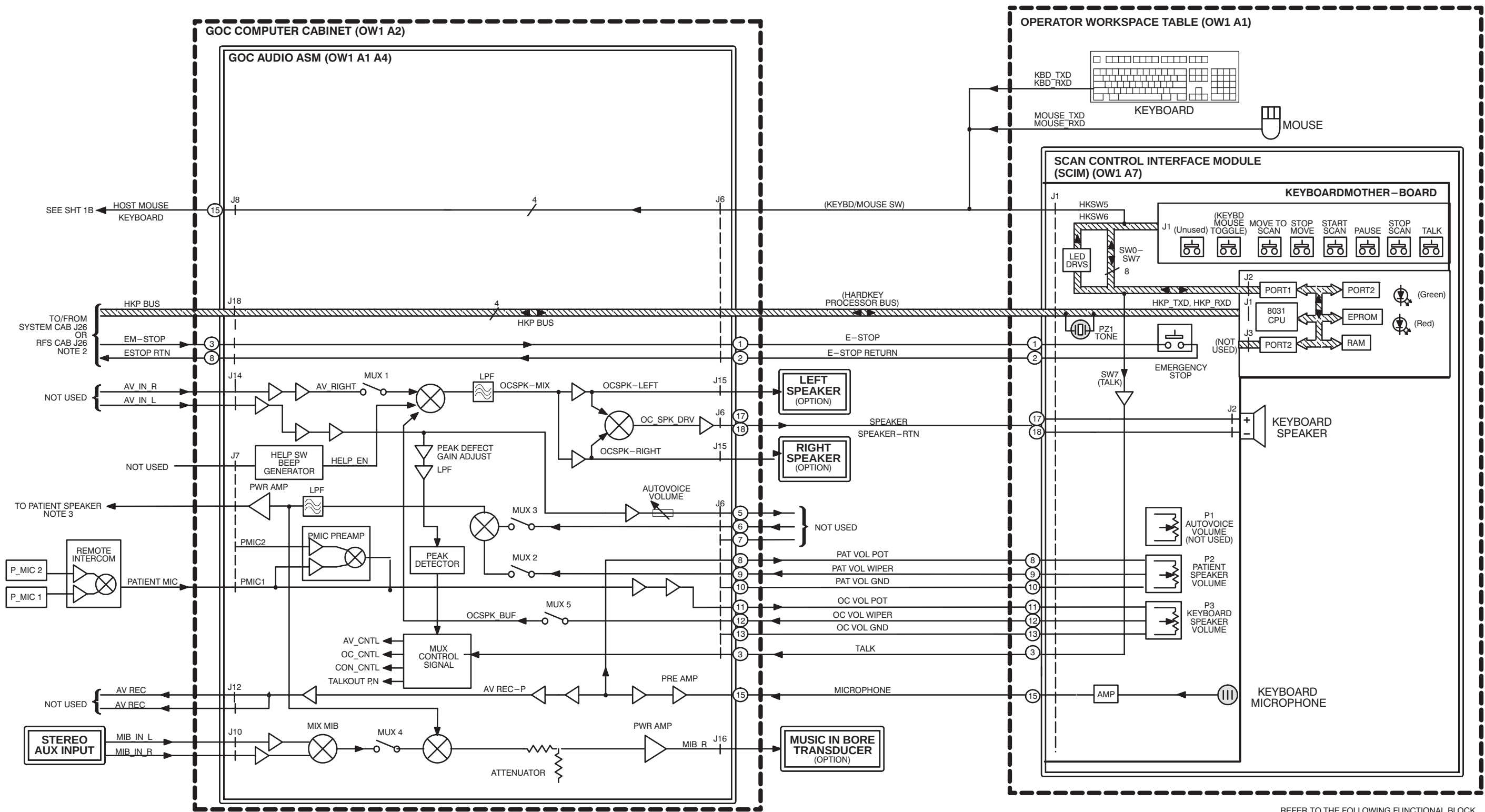


G3 & 3T94 EXCITE HD (12.x) OVERALL SYSTEM
(OPERATOR WORKSPACE CABINET/TABLE with GOC AUDIO ASM)

SHEET 2B OF 2

REFER TO THE FOLLOWING FUNCTIONAL BLOCK
DIAGRAM FOR ADDITIONAL INFORMATION:

NOTE:
1) PATIENT HANDLING
2) MGD/RRF/RF
3) OVERALL SYSTEM

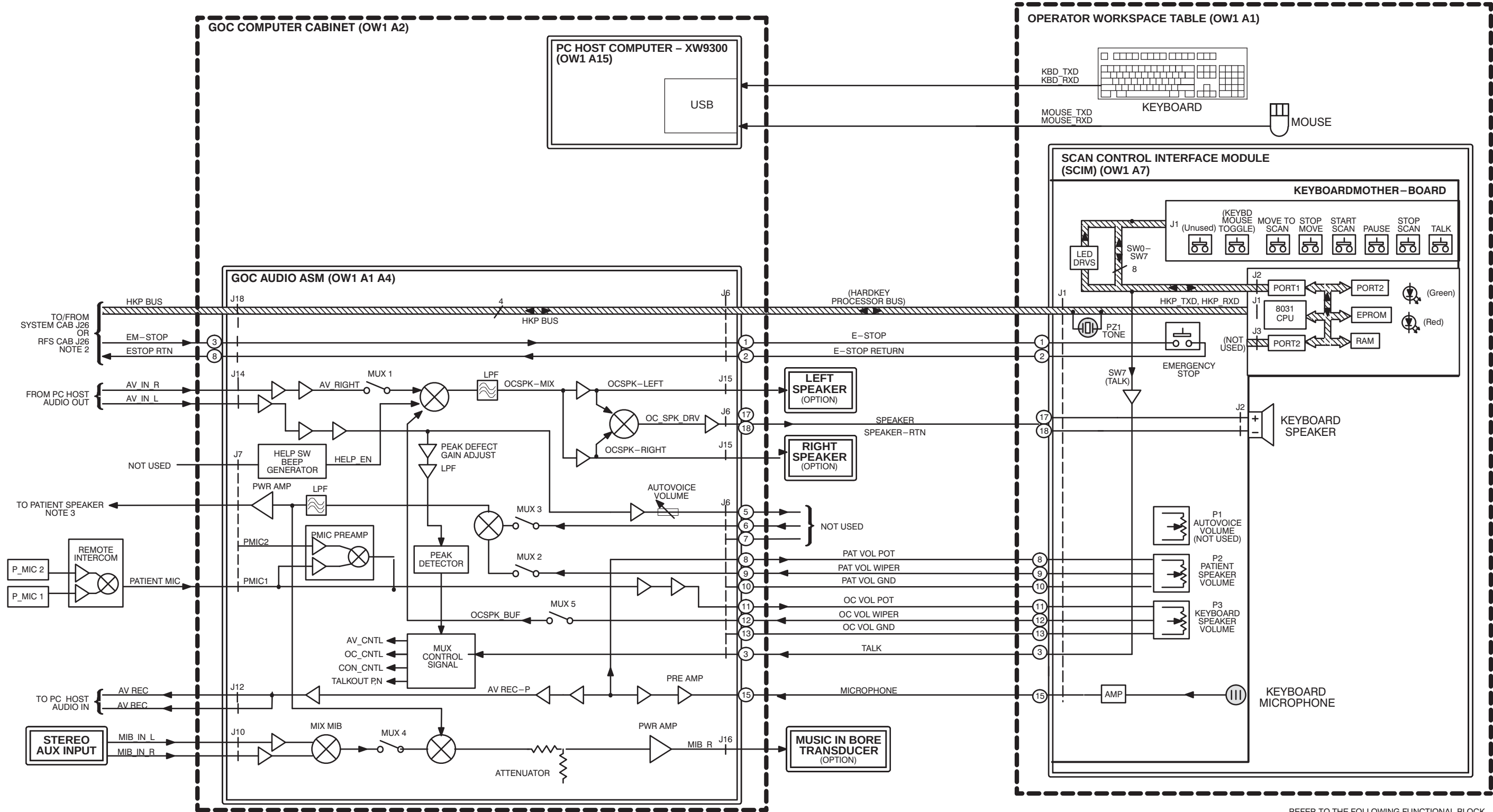

G3 & 3T94 HDx (14.x) & HDxT (15.x) OVERALL SYSTEM
(OPERATOR WORKSPACE CABINET/TABLE with GOC AUDIO ASM)

SHEET 2C OF 2


REFER TO THE FOLLOWING FUNCTIONAL BLOCK
DIAGRAM FOR ADDITIONAL INFORMATION:


NOTE:
1) PATIENT HANDLING
2) MGD/RRF/RF
3) OVERALL SYSTEM


PATIENT HANDLING




TABLE OF CONTENTS

<u>SECTION</u>	<u>PAGE</u>
1 – SYSTEM	
G3 & 3T94 HD, HDx & HDxt PATIENT HANDLING (Forward Production or Pre-EXCITE Upgrades)	1-1
G3 & 3T94 HD & HDx & HDxt PATIENT HANDLING (EXCITE Upgrades)	1-2
G3 & 3T94 HD & HDx & HDxt PATIENT HANDLING (Forward Production or Pre-EXCITE Upgrades)	1-3
G3 & 3T94 HD & HDx & HDxt PATIENT HANDLING (EXCITE Upgrades)	1-4
G3 & 3T94 HD & HDx & HDxt PATIENT HANDLING POWER SUPPLY (PHPS)	1-5
G3 HD & HDx & HDxt PATIENT HANDLING – SSM (EXCITE Upgrades)	1-6
G3 HD & HDx & HDxt PATIENT HANDLING – SSM (EXCITE Upgrades)	1-7
3T94 HD & HDx & HDxt PATIENT HANDLING – SSM (EXCITE Upgrades)	1-8
3T94 HD & HDx & HDxt PATIENT HANDLING – SSM (EXCITE Upgrades)	1-9

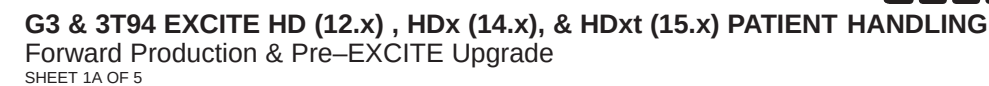


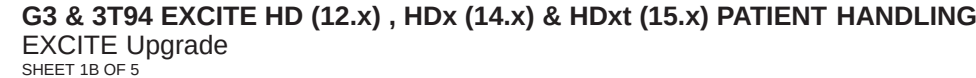


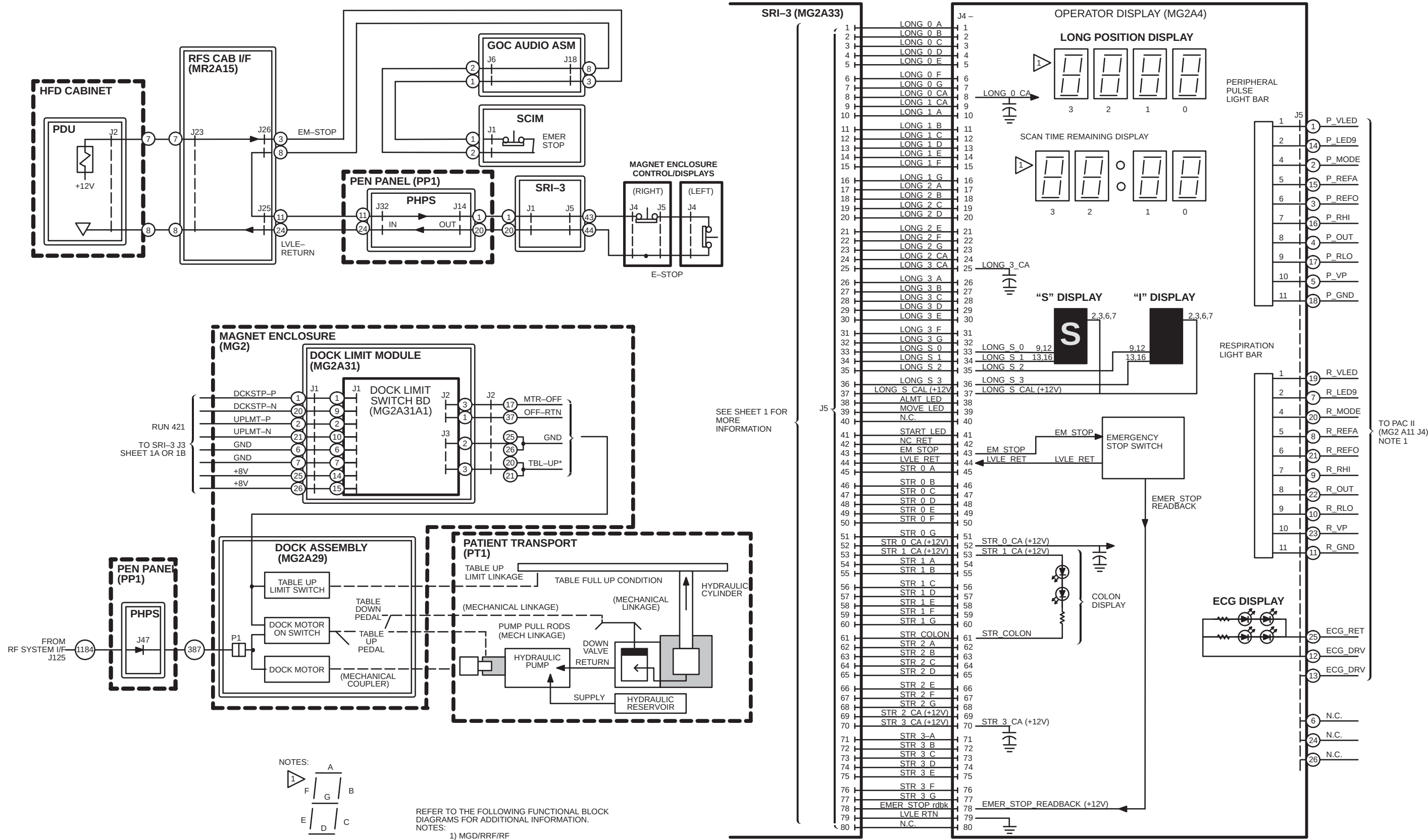

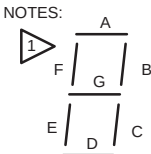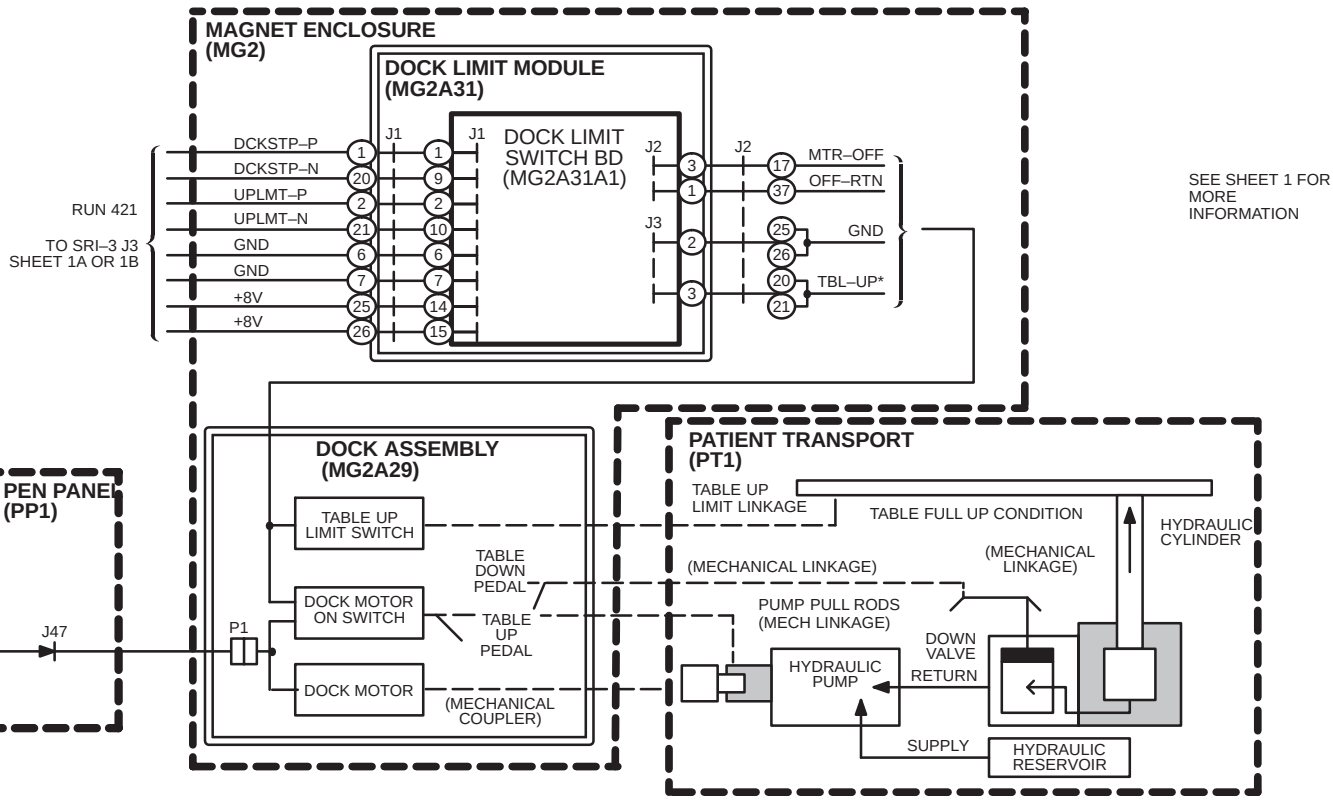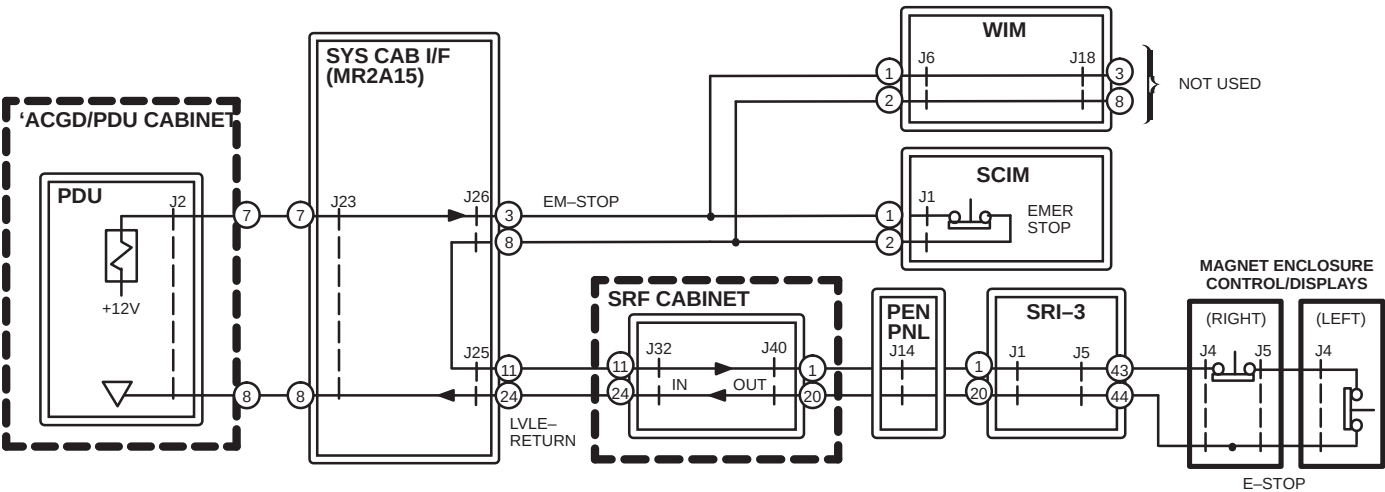


REFER TO THE FOLLOWING FUNCTIONAL BLOCK
DIAGRAMS FOR ADDITIONAL INFORMATION.
NOTES:
1) MGD/RRF/RF

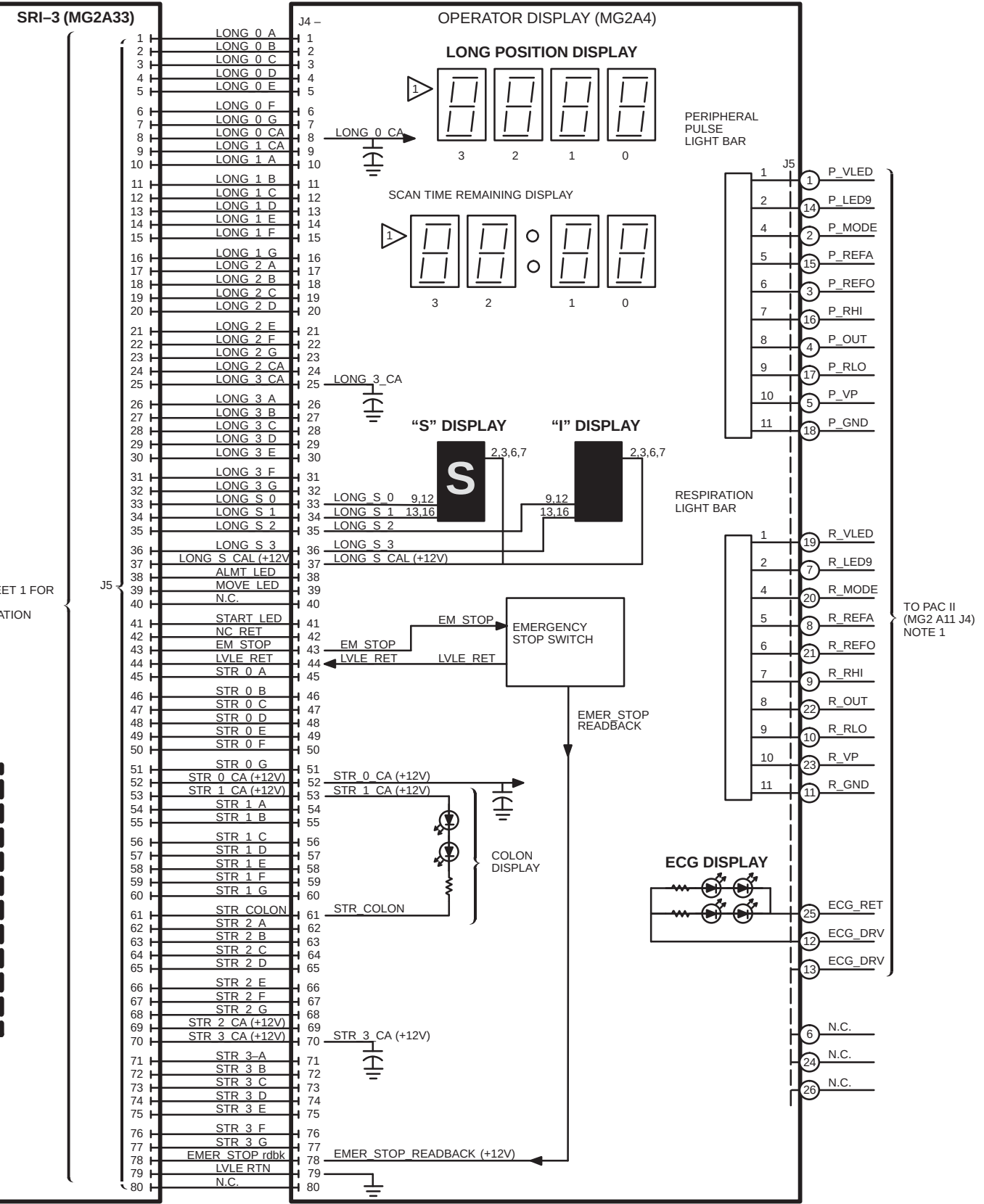

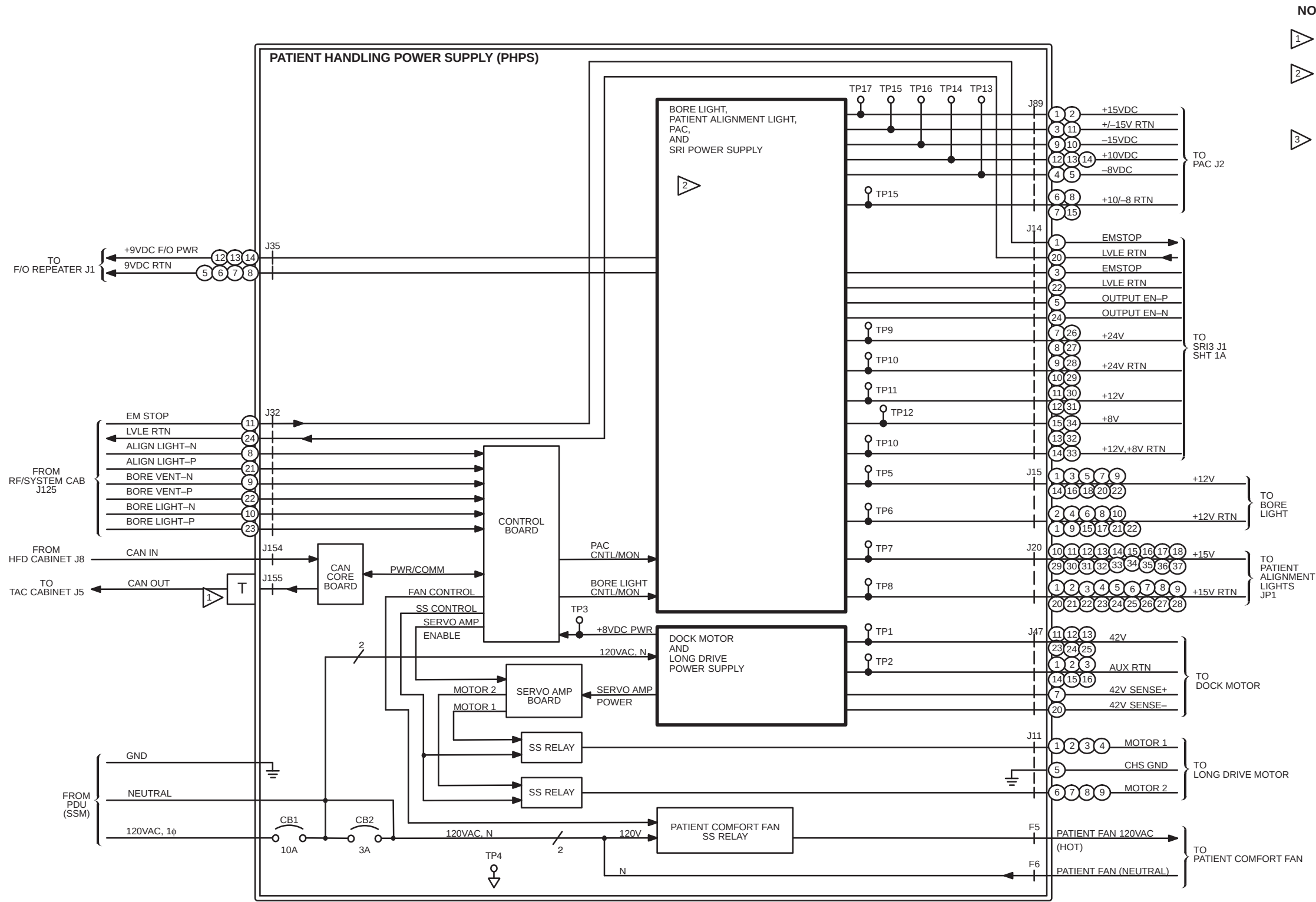


- NOTES:**
- 1 OPTIONAL: Used only if High Order Shim (HOS) installed. PP1 PHPS J155 Terminator removed when HOS is installed.
 - 2 Systems equipped with Laser Patient Alignment lights use one type of power supply while those equipped with Halogen patient Alignment lights must use another. See the Renewal Parts manual.
 - 3

TEST POINT	FUNCTION
TP1	TABLE VOLTS (42VDC)
TP2	TABLE COM
TP3	CONTROL BD +8V
TP4	+8V COM
TP5	BORE +12V
TP6	BORE COM
TP7	PAL +15V
TP8	PAL COM
TP9	SRI +24V
TP10	SRI COM
TP11	SRI +12V
TP12	SRI +8V
TP13	PAC -8V
TP14	PAC +10V
TP15	PAC COM
TP16	PAC -15V
TP17	PAC +15V

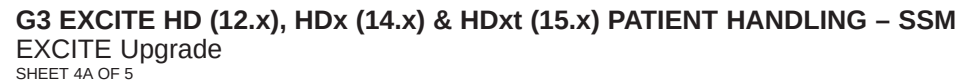


REFER TO THE FOLLOWING FUNCTIONAL BLOCK
DIAGRAM FOR ADDITIONAL INFORMATION:
NOTE:
1) MGD/RRF/RF



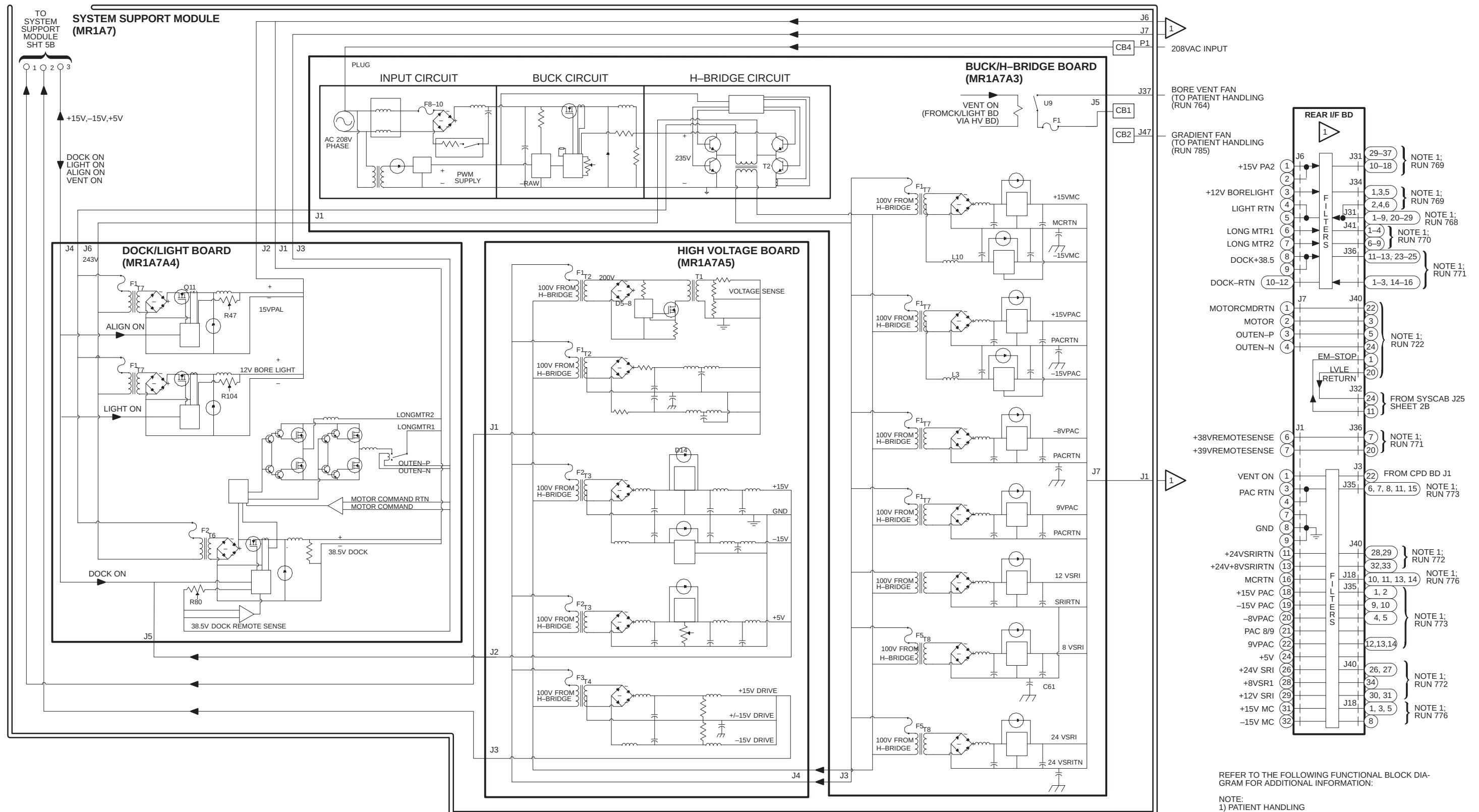


3T94 EXCITE HD (12.x), HDx (14.x) & HDxt (15.x) PATIENT HANDLING – SSM
EXCITE Upgrade
SHEET 5A OF 5

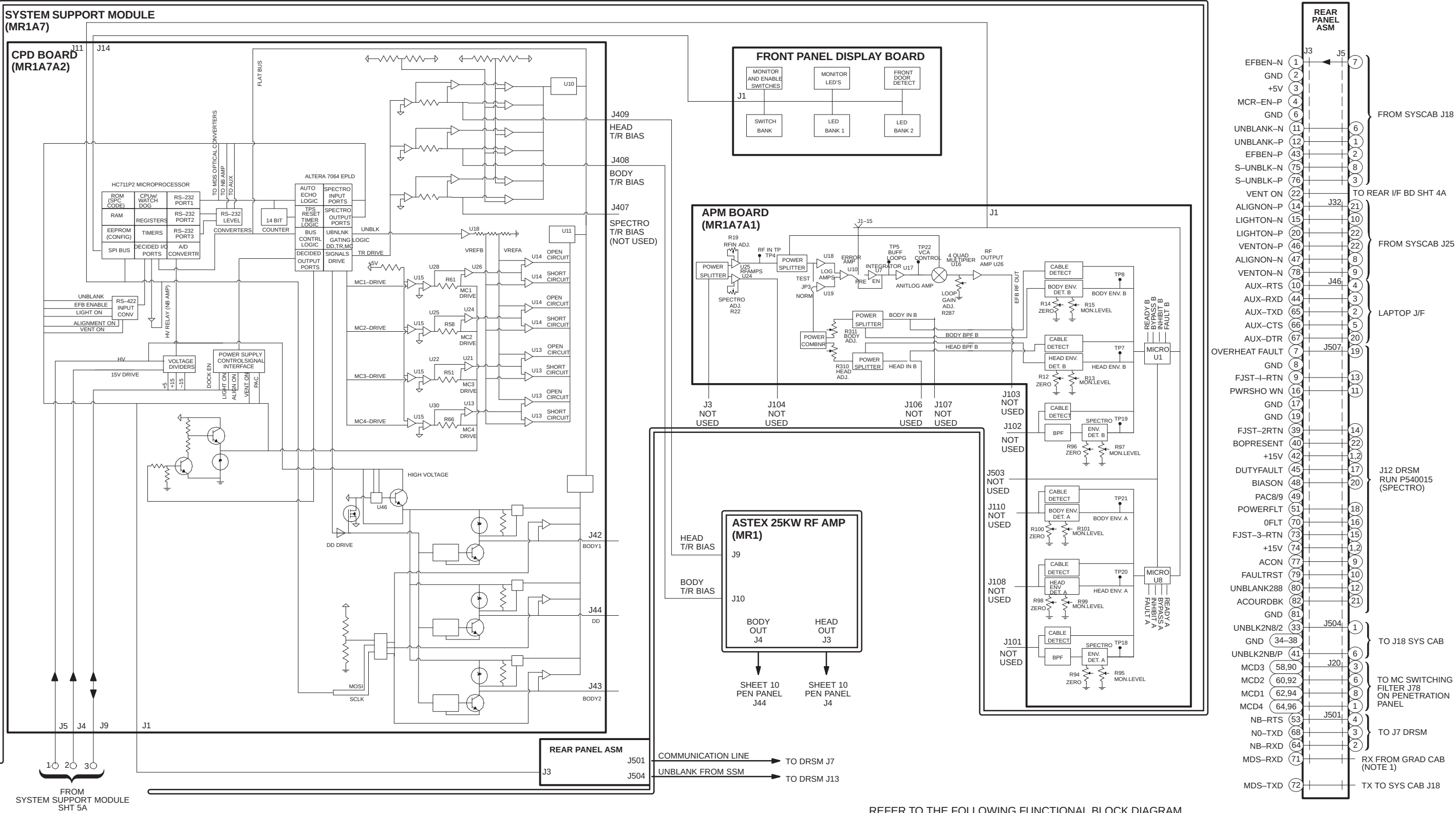

Refer to the following Functional Block Diagrams for additional information:

1 OPTIONAL: Used only if Patient Handling Power Supply (PHPS) is not installed. J8 terminator removed when PHPS is installed.

2 Present only on Forward Production or Pre-EXCITE upgrade systems.

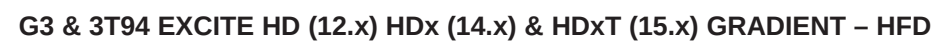


NOTE:
Additional Signal/Cabling detail can be found under Interconnects on the Service Methods CDROM


Refer to the following Functional Block Diagrams for additional information:

- 1) MGD/RRF/RF
- 2) PATIENT HANDLING

- 1 Terminator present at J8 if system is not a TwinSpeed.
- 2 TAC Cabinet present only on TwinSpeed systems.

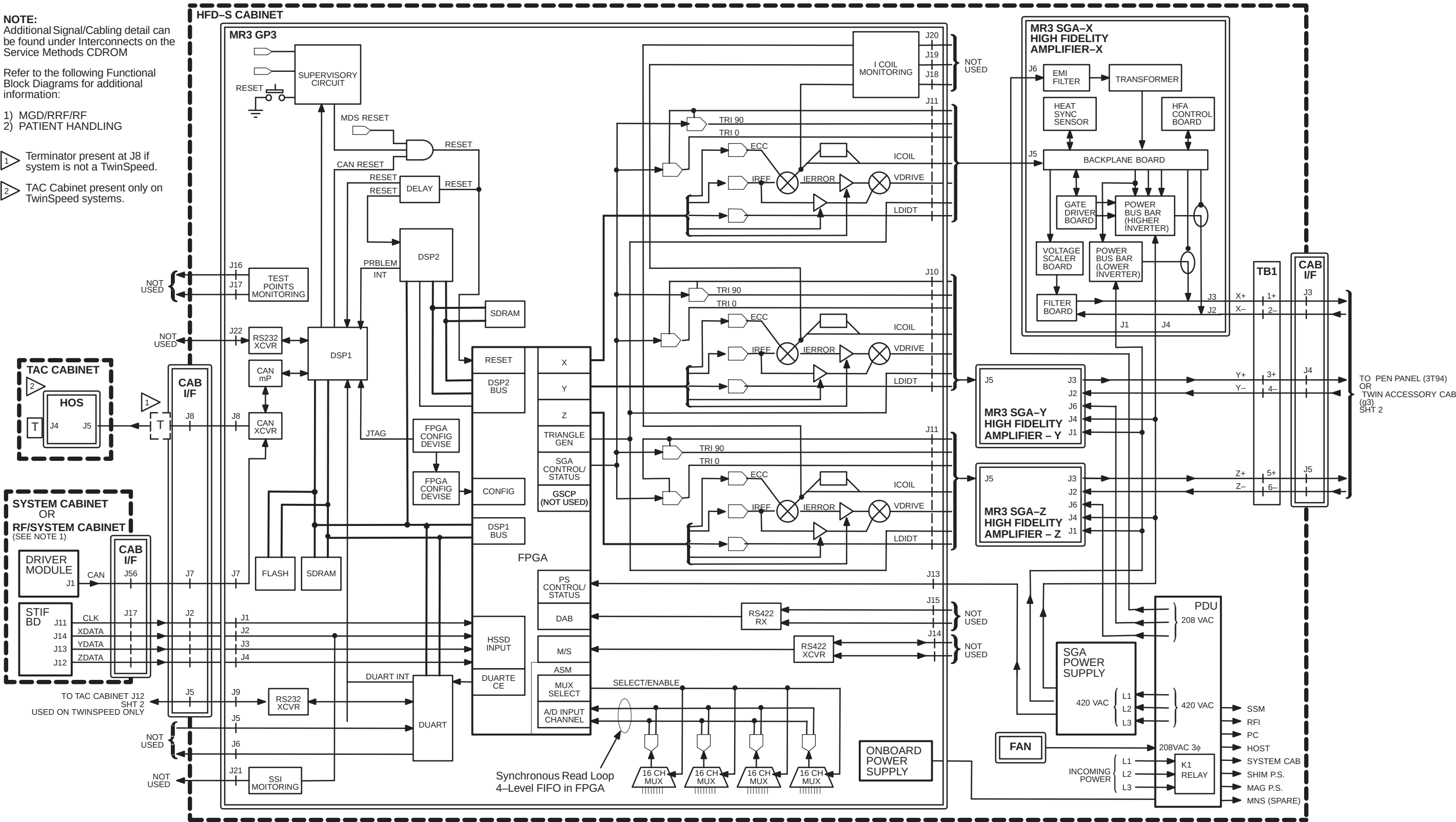

G3 & 3T94 EXCITE HD (12.x), HDx (14.x) & HDxt (15.x) GRADIENT – HFD-S (Upgraded ACGD)

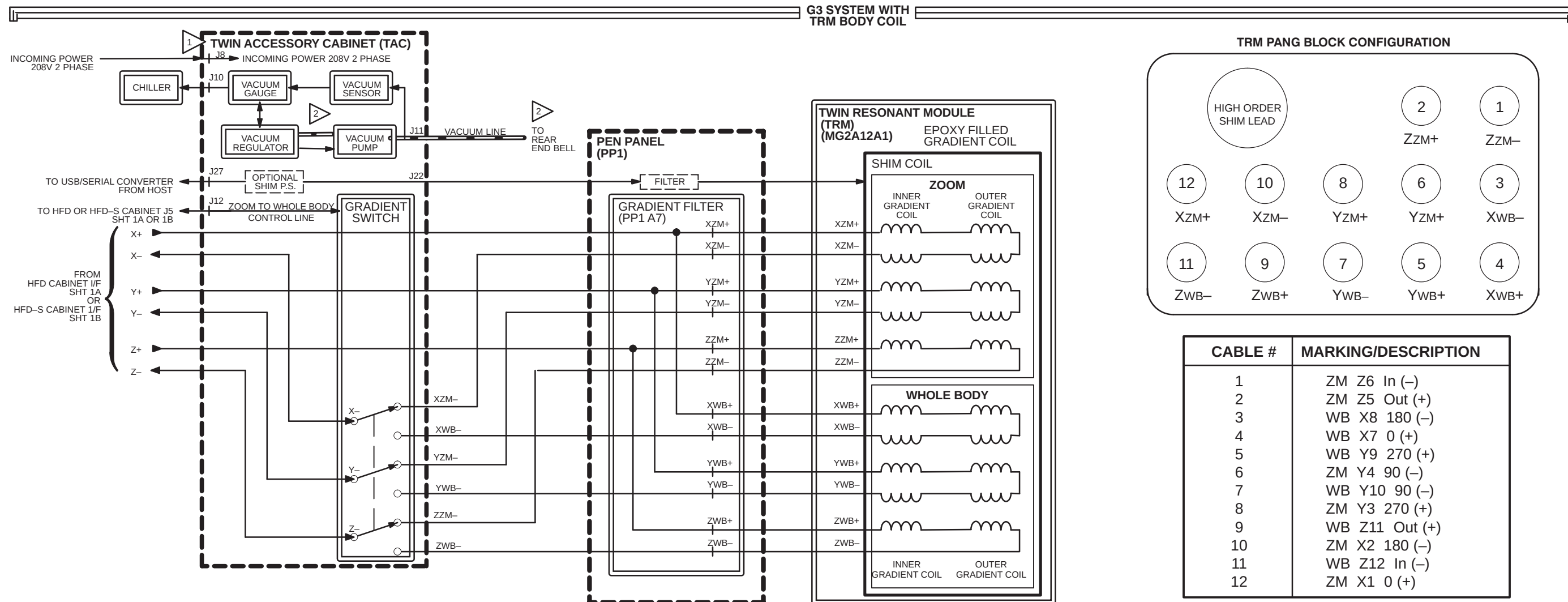


NOTE:
Refer to the following Functional Block Diagrams for additional information:

1) OVERALL SYSTEM

- 1 Twin Accessory Cabinet is present only on TwinSpeed systems (with TRM body coil).
- 2 Vacuum system is not installed for HDx (14.x) systems. It is used only on HD (12.x) systems.

GRADIENT CHAIN HARDWARE

SHEET 2 OF 2

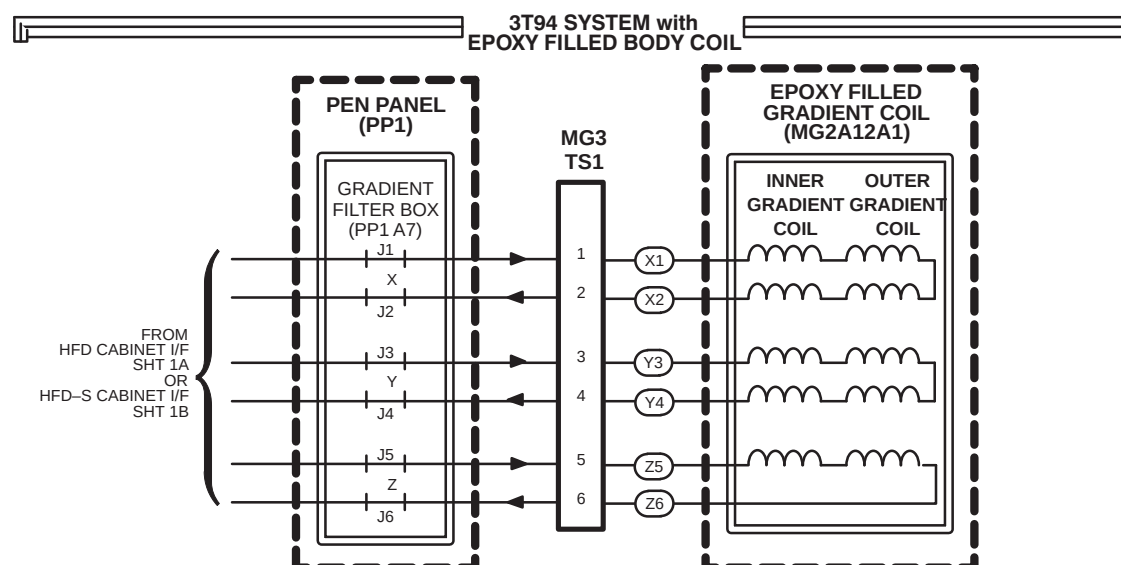

PDU

TABLE OF CONTENTS

Note:

PDU found in Vendor Manuals only. (SEE SERVICE METHODS CDROM)

MGD/RRF/RF

TABLE OF CONTENTS


<u>SECTION</u>	<u>PAGE</u>
1 – MGD/RRF/RF SUBSYSTEM	
G3 & 3T94 HD MGD CHASSIS – HOST, RECEIVE, AND RECON FUNCTIONS (Forward Production & Upgrades)	1–1
G3 & 3T94 HDx & HDxt MGD CHASSIS – HOST, RECEIVE, AND RECON FUNCTIONS (Forward Production & Upgrades)	1–2
G3 & 3T94 HDx & HDxt VOLUME RECON ENGINE (Forward Production & Upgrades)	1–3
G3 & 3T94 HD, HDx & HDxt MGD CHASSIS – SEQUENCE SEGMENT	1–4
G3 & 3T94 HD, HDx & HDxt RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 1 OF 2A (8 Channel System)	1–5
G3 & 3T94 HD, HDx & HDxt RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 1 OF 2B (16 Channel – Forward Production & Upgrades)	1–6
G3 HDx & HDxt RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 1 OF 4C (32 Channel – Forward Production & Upgrades)	1–7
G3 & 3T94 HD, HDx & HDxt RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 2 OF 2A (8 Channel System)	1–8
G3 & 3T94 HD, HDx & HDxt RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 2a OF 2B (Channels 1–8 of 16 Channel – Forward Production & Upgrades)	1–9
G3 & 3T94 HD, HDx & HDxt RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 2a OF 4C (Channels 1–8 of 32 Channel – Forward Production & Upgrades)	1–10
G3 & 3T94 HD, HDx & HDxt RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 2b OF 2B (Channels 9–16 of 16 or 32 Channel – Forward Production & Upgrades)	1–11
G3 HDx & HDxt RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 2c OF 4C (Channels 17–24 of 32 Channel – Forward Production & Upgrades)	1–12
G3 HDx & HDxt RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 2d OF 4C (Channels 25–32 of 32 Channel – Forward Production & Upgrades)	1–13
G3 & 3T94 HD, HDx & HDxt RRF RECEIVE CHAIN CONNECTIONS	1–14
G3 HD DRIVER MODULE (Upgrade with SSM)	1–15
G3 HD DRIVER MODULE (Forward Production, RFS Cabinet, no SSM)	1–16
G3 HDx & HDxt DRIVER MODULE (Upgrade with SSM)	1–17
G3 HDx & HDxt DRIVER MODULE (Forward Production, RFS Cabinet, no SSM)	1–18
3T94 HDx & HDxt DRIVER MODULE (Upgrade with SSM)	1–19
3T94 HDx & HDxt DRIVER MODULE (Pre-EXCITE Upgrade, RFS Cabinet, no SSM)	1–20
G3 HD, HDx & HDxt SYSTEM SUPPORT MODULE (Upgrades Only)	1–21
3T94 HD, HDx & HDxt SYSTEM SUPPORT MODULE (Upgrades Only)	1–22


TABLE OF CONTENTS (Continued)

<u>SECTION</u>	<u>PAGE</u>
1 – MGD/RRF/RF SUBSYSTEM	
G3 HD and HDx & HDxt 35kW MKS RF AMPLIFIER	1–23
G3 & 3T94 HD, HDx & HDxt MAGNET ROOM & RELATED COMPONENTS (HD – Forward Production & Upgrades, HDx – Upgrades)	1–24
G3 HDx & HDxt MAGNET ROOM & RELATED COMPONENTS (Forward Production)	1–25
3T94 HD, HDx & HDxt MAGNET ROOM & RELATED COMPONENTS (Upgrades)	1–26
G3 & 3T94 HD, HDx & HDxt MULTI-COIL MODULES (HD – Forward Production & Upgrades, HDx – Upgrades)	1–27
G3 HDx & HDxt 8 and 16 CHANNEL MULTI-COIL	1–28
G3 HDx & HDxt 32 CHANNEL MULTI-COIL	1–29
G3 HDx & HDxt LPCA MULTI-COIL INTERFACE ASM	1–30
G3 & 3T94 HD, HDx & HDxt HYPERTRONICS CONNECTOR PINOUTS (Rows A & B) (HD – Forward Production & Upgrades, HDx – Upgrades)	1–31
G3 & 3T94 HD, HDx & HDxt HYPERTRONICS CONNECTOR PINOUTS (Rows C, D, & E) (HD – Forward Production & Upgrades, HDx – Upgrades)	1–32
G3 HD and HDx & HDxt RF ACCESSORIES CABINET (Upgrades)	1–33
3T94 HD RF ACCESSORIES CABINET (Upgrades)	1–34
G3 & 3T94 HD, HDx & HDxt PHYSIOLOGICAL ACQUISITION CONTROLLER II	1–35
G3 & 3T94 HD, HDx & HDxt 16 CHANNEL SWITCH (3.0T)	1–36
HDx & HDxt DEPOPULATED 8 CHANNEL SWITCH (3.0T)	1–37
G3 & 3T94 HD 8 CHANNEL MULITCOIL “A” AND “B” CONNECTOR	1–38
G3 HDx & HDxt MULITCOIL “A+” CONNECTOR (Forward Production)	1–39
G3 HDx & HDxt MULITCOIL “B8” CONNECTOR (Forward Production)	1–40
G3 HDx & HDxt MULITCOIL “B16” and “C16” CONNECTOR (Forward Production)	1–41
G3 & 3T94 HD, HDx & HDxt ASC CHASSIS	1–42
G3 HD, HDx & HDxt MNS INTERCONNECTS	1–43
G3 & 3T94 HD, HDx & HDxt MCR II TOOL AND CABLE (HD – Forward Production & Upgrades, HDx – Upgrades)	1–44
G3 HDx & HDxt MCR III TOOL AND CABLE (Forward Production)	1–45
G3 & 3T94 HD, HDx & HDxt MCR II TOOL – 8 CHANNEL CONNECTOR (HD – Forward Production & Upgrades, HDx – Upgrades)	1–46
G3 HDx & HDxt MCR III TOOL – 8 CHANNEL CONNECTOR (Forward Production)	1–47
2 – MGD/RRF/RF SUBSYSTEM SUPPLEMENTAL SCHEMATICS	
DYNAMIC DISABLE FILTER MODULE SCHEMATIC	2–1

NOTE:

Additional Signal/Cabling detail can be found under Interconnects on the Service Methods CDROM.





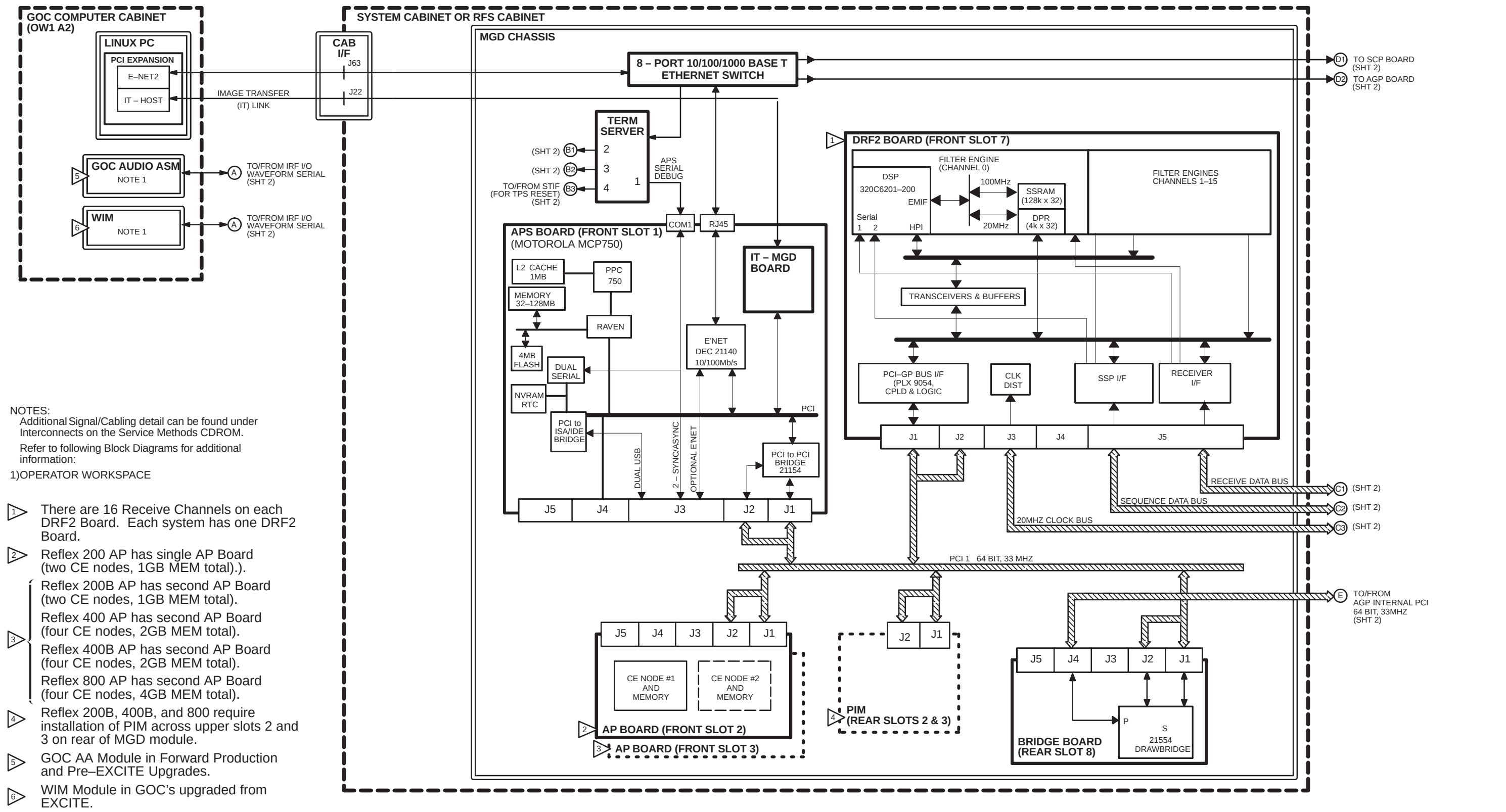


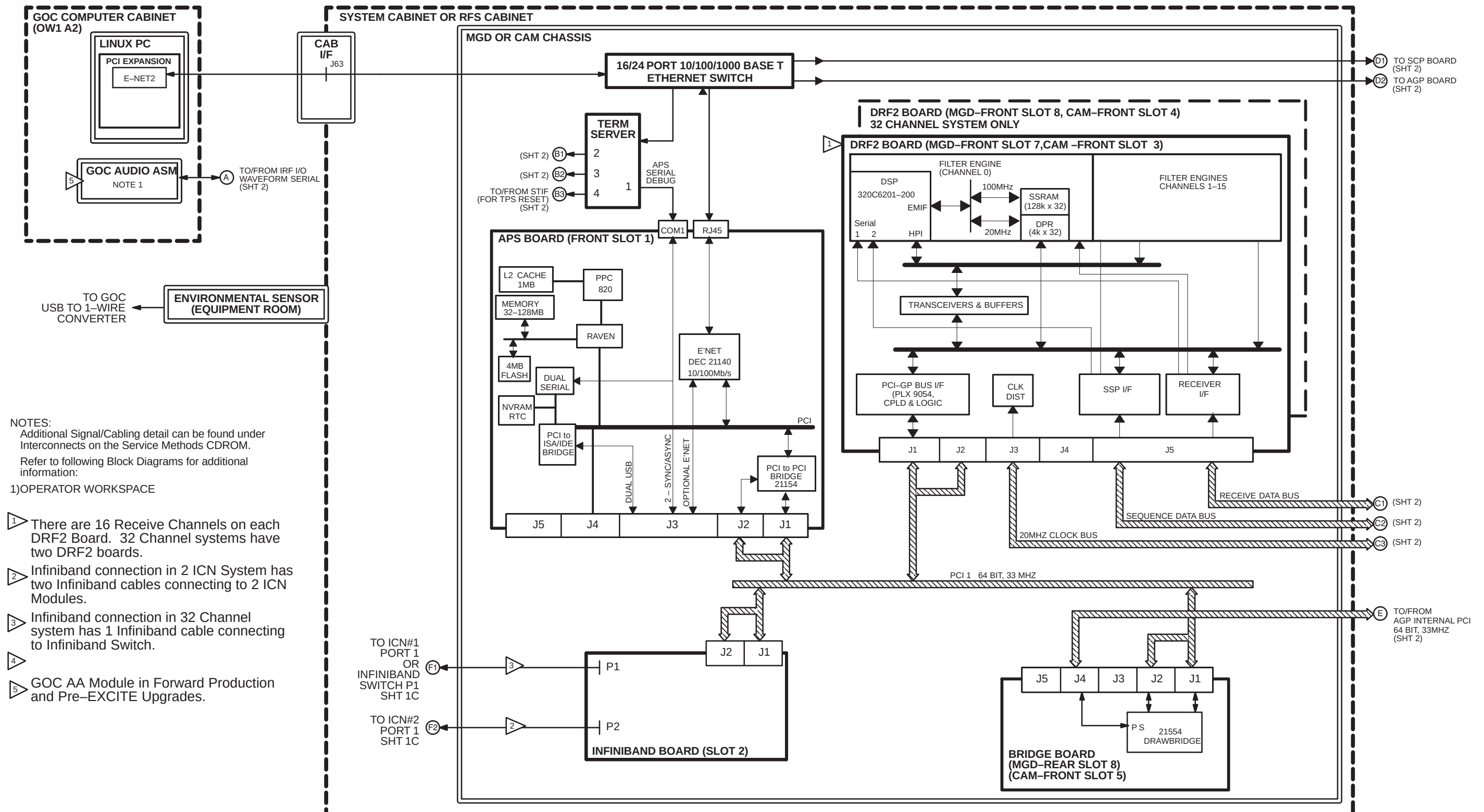



NOTES:
Additional Signal/Cabling detail can be found under Interconnects on the Service Methods CDROM.
Refer to following Block Diagrams for additional information:

1) OPERATOR WORKSPACE

Other connections to 16/24 port Ethernet Switch on sheet 1B.

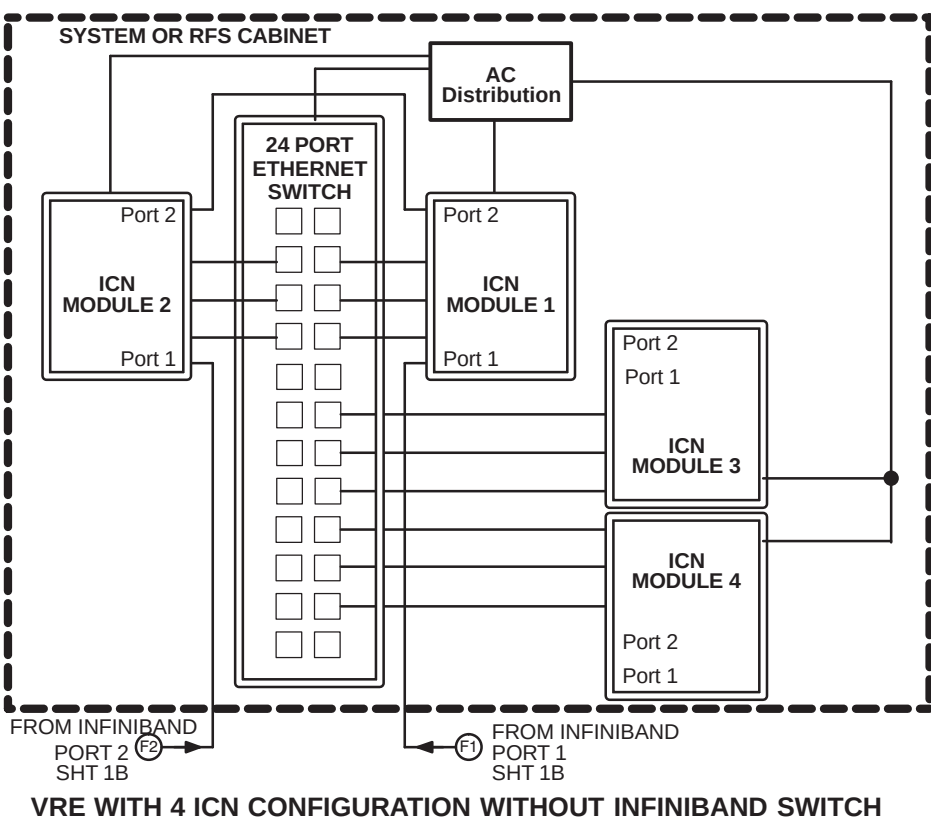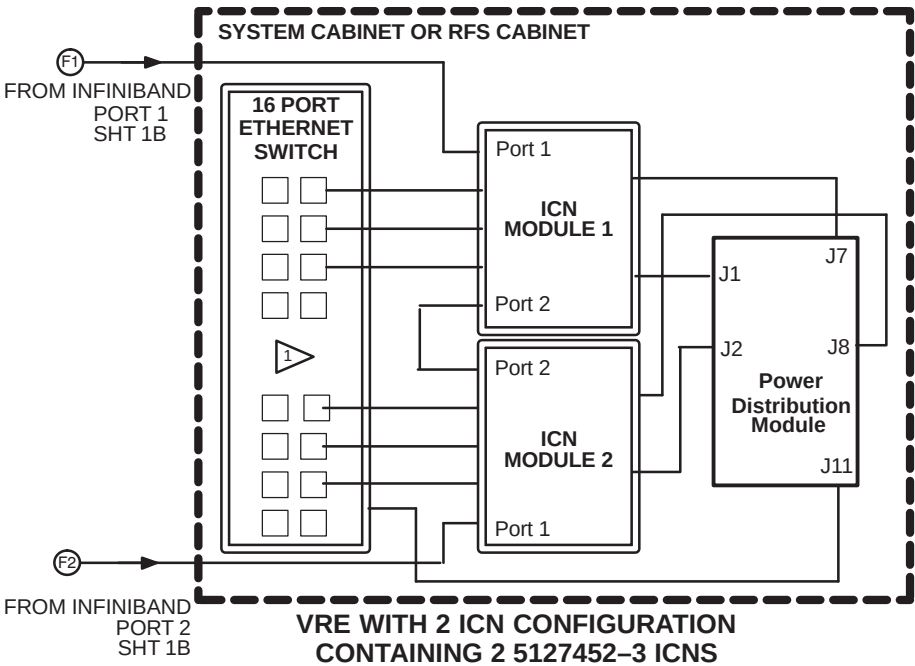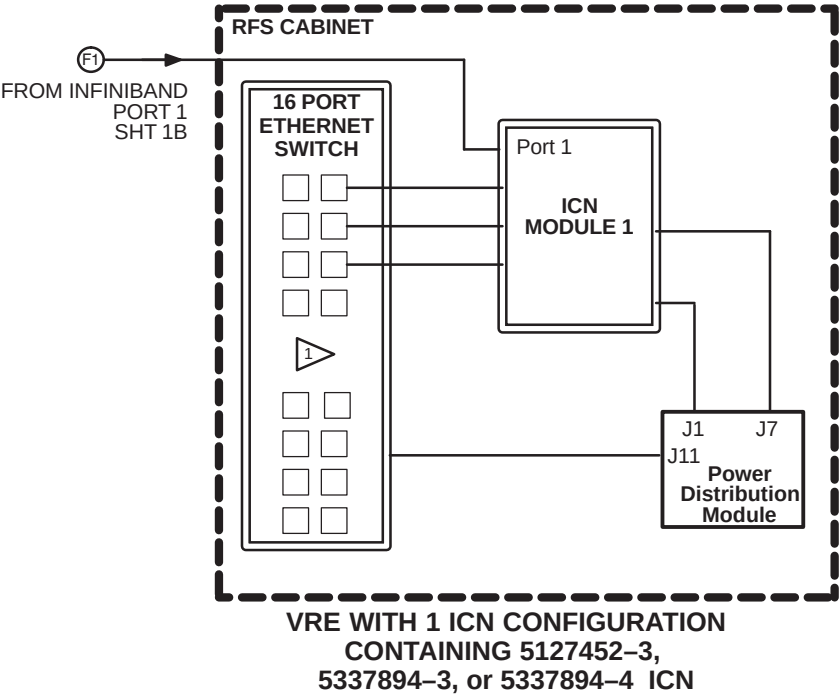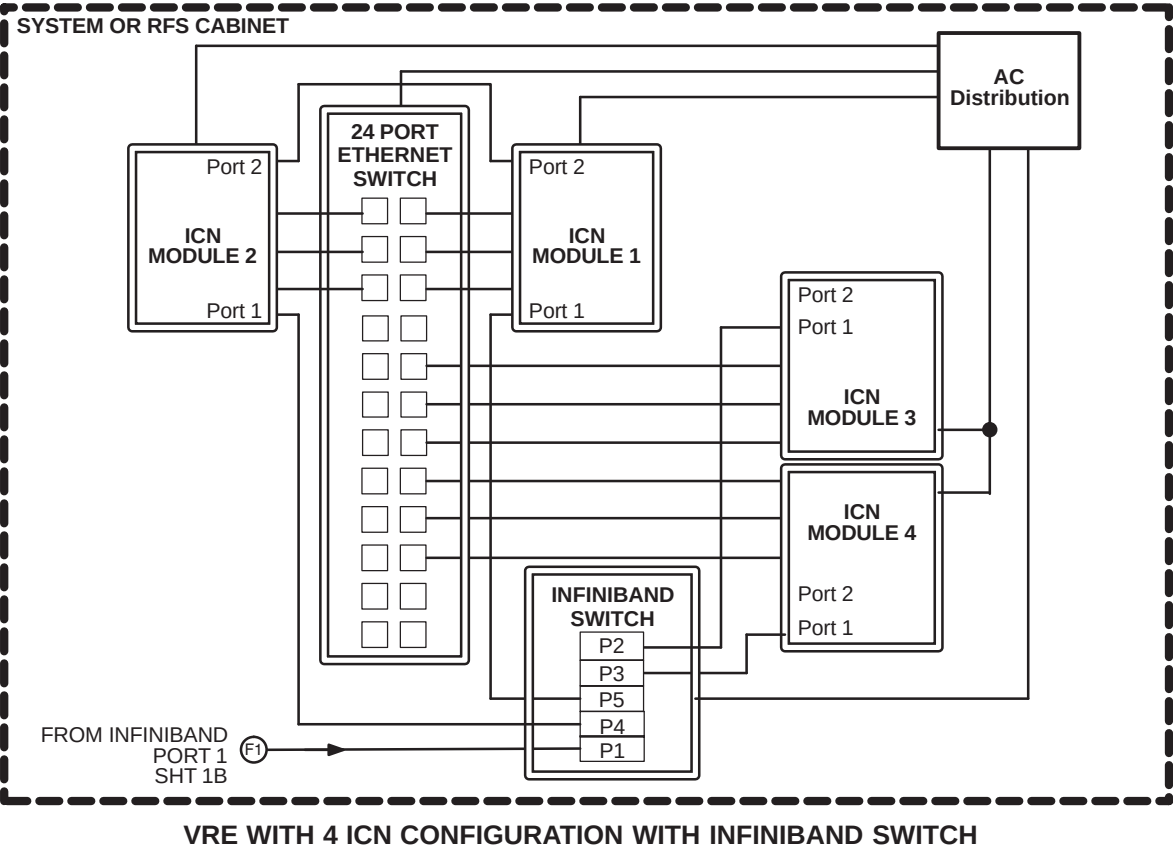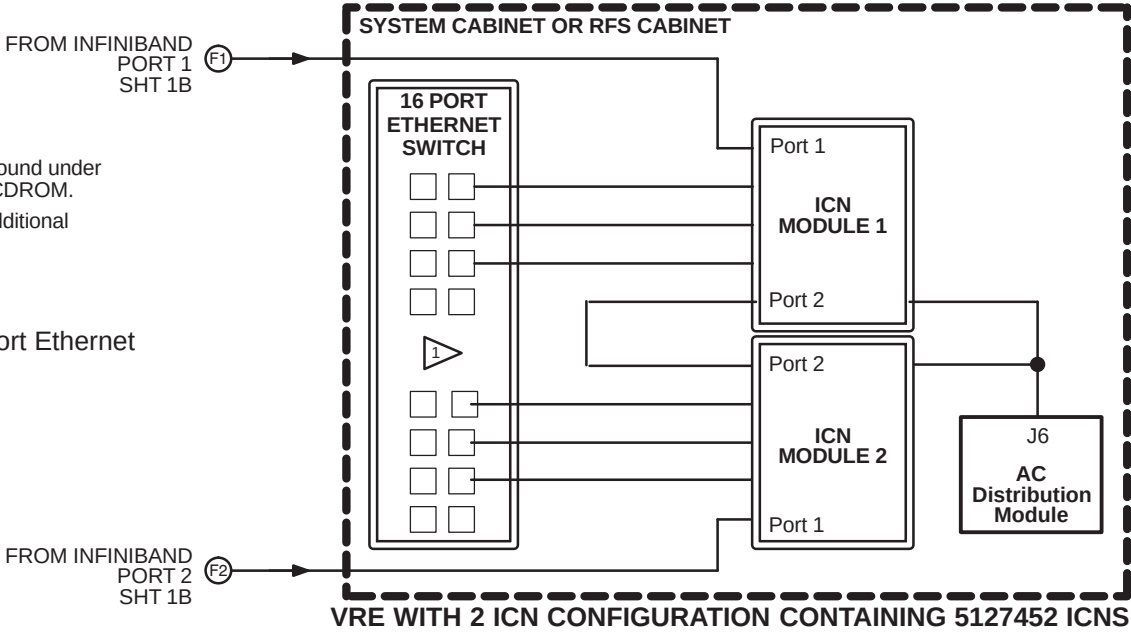


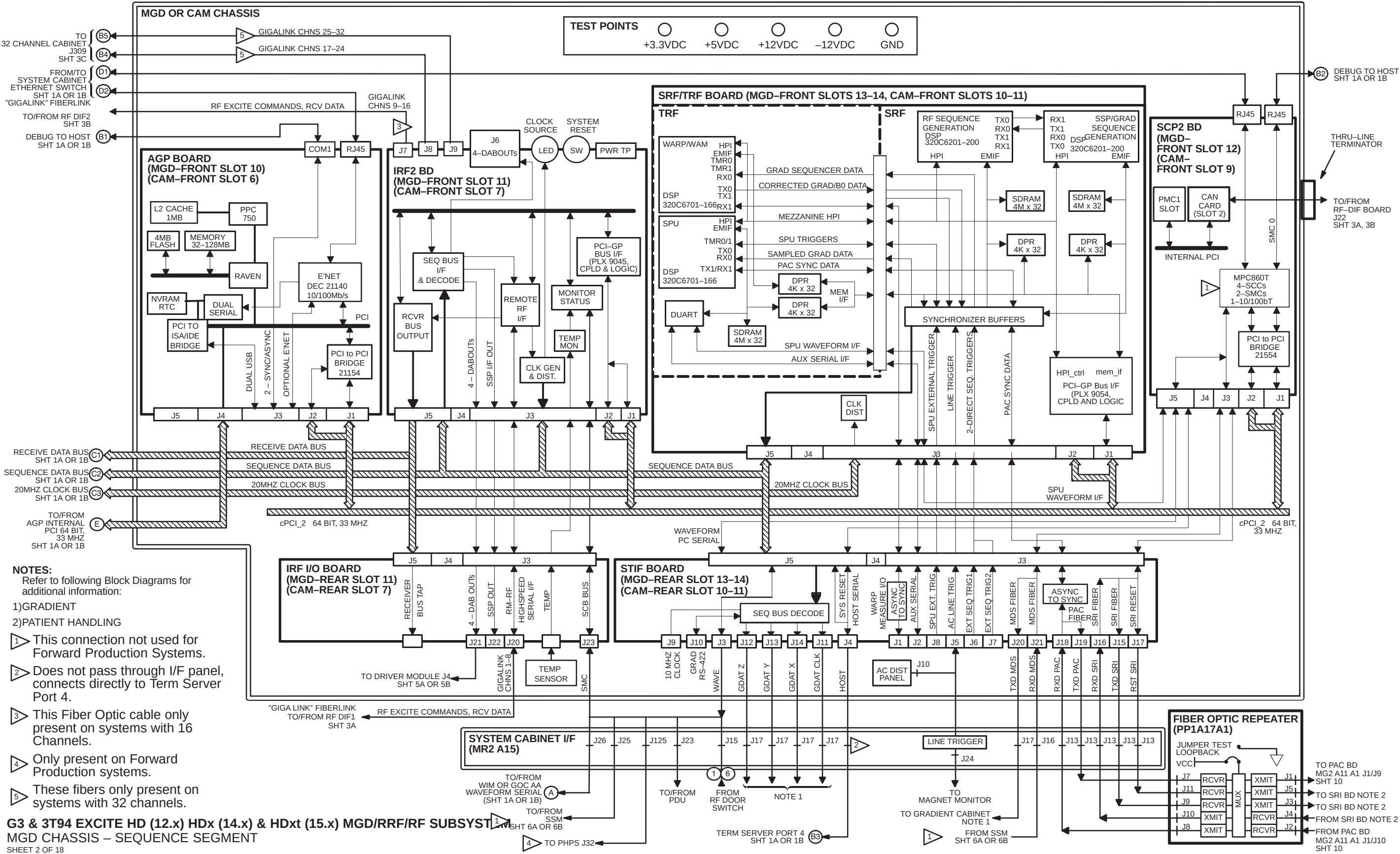

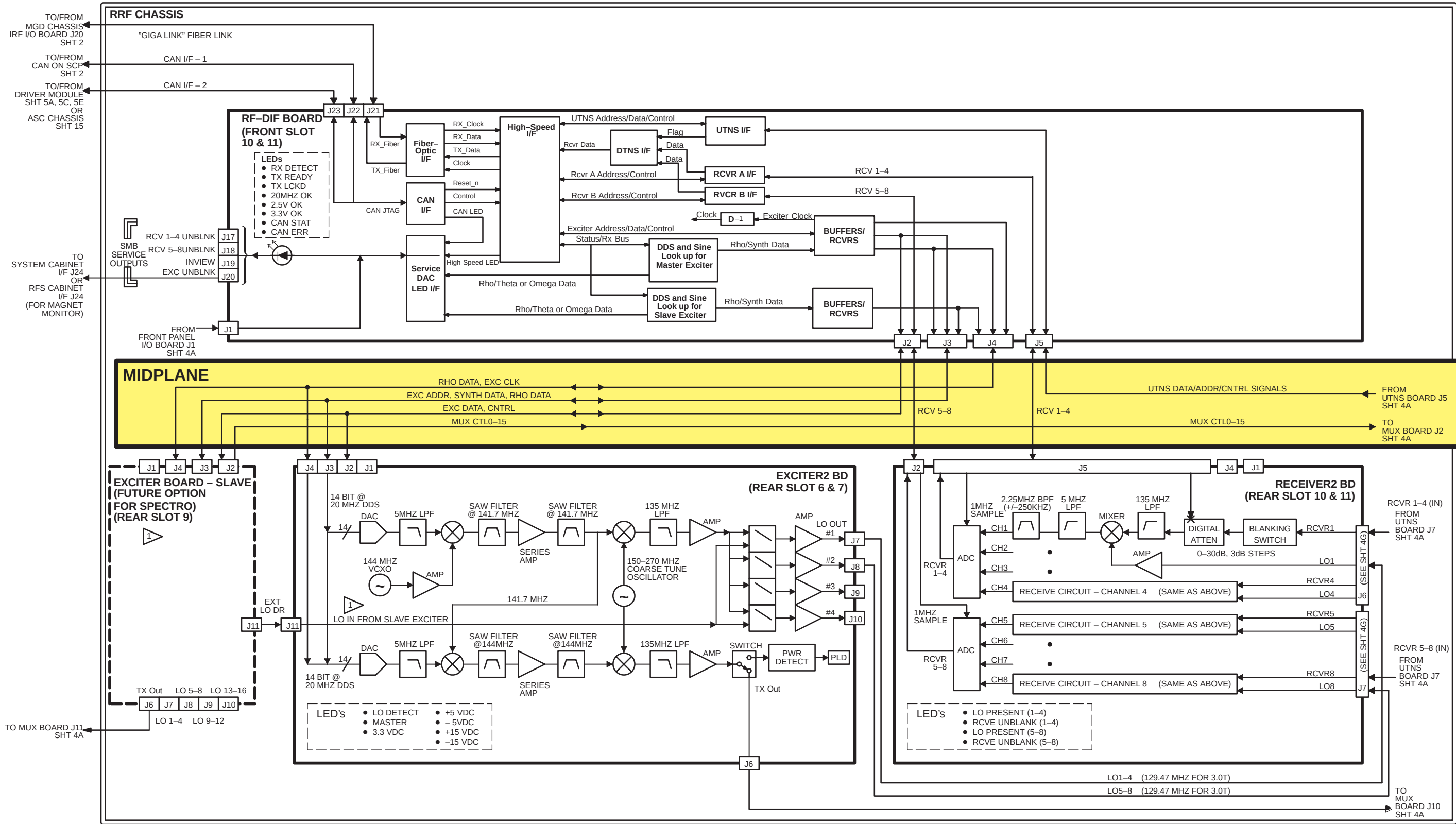

G3 & 3T94 EXCITE HD (12.x), HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
(8 Channel Systems)
RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 1 OF 2A (8 Channel)
SHEET 3A OF 18

NOTE:
➤ Slave Exciter is a future option for Spectro.

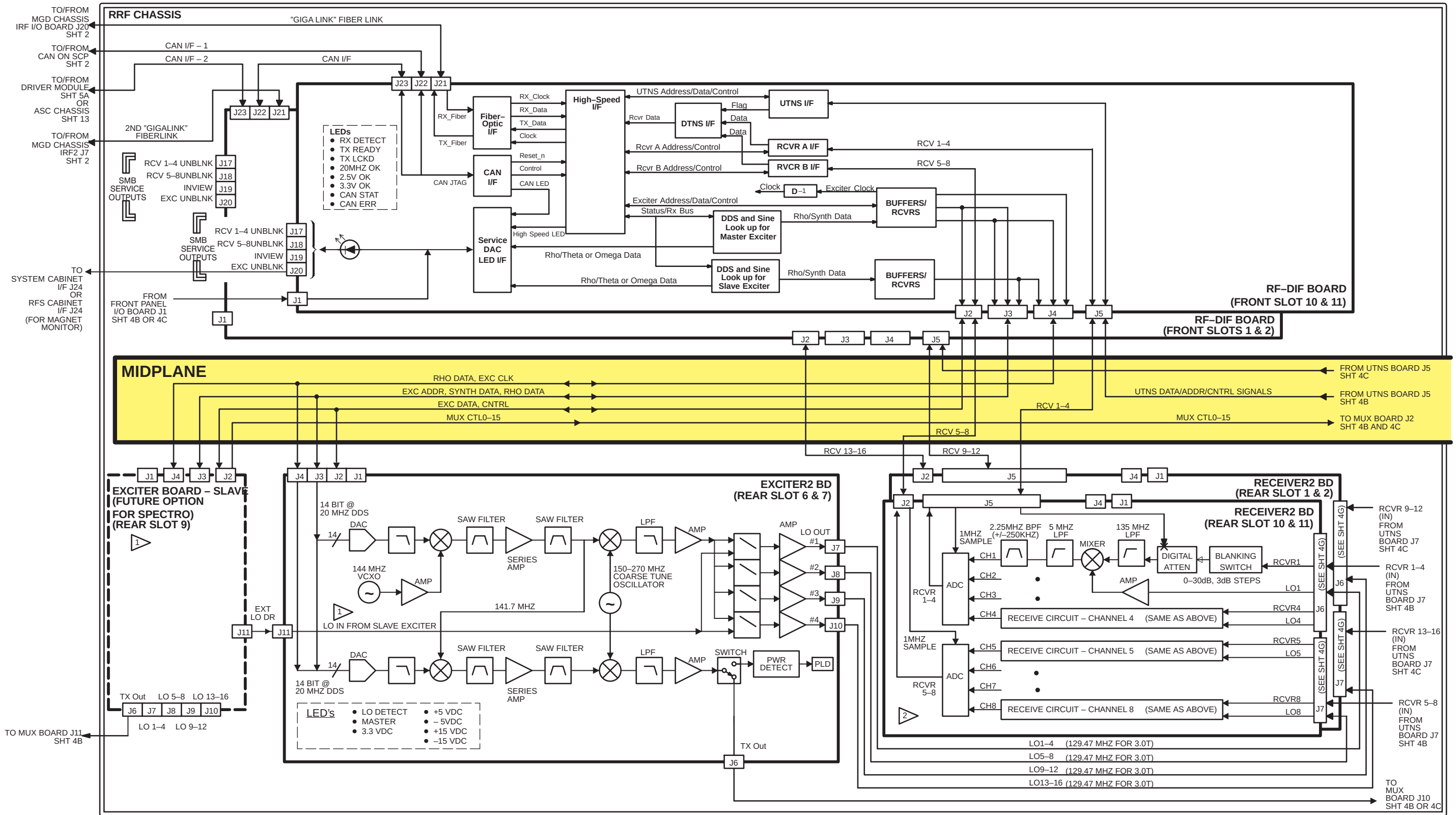


G3 & 3T94 EXCITE HD (12.x) , HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
(16 Channel – Forward Production & Upgrades)
RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 1 OF 2B (16 Channel)
SHEET 3B OF 18

NOTE:
1 Slave Exciter is a future option for Spectro.

2 LED's


- LO PRESENT (1-4)
- RCVE UNBLANK (1-4)
- LO PRESENT (5-8)
- RCVE UNBLANK (5-8)

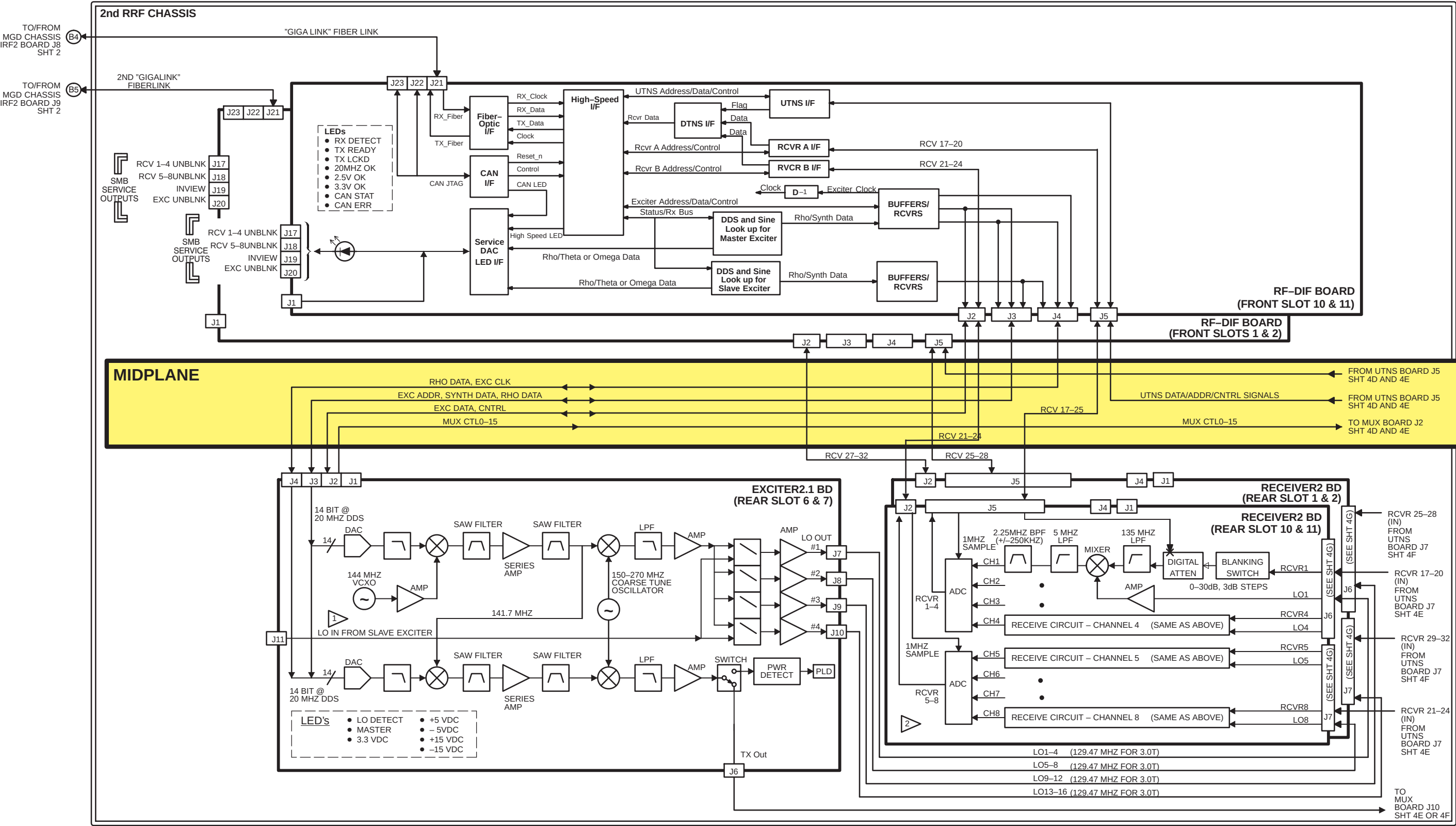

G3 HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
(32 Channel – Forward Production & Upgrades)
2nd RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 1 OF 4C (32 Channel)
SHEET 3C OF 18


NOTE: **LED's**

- LO PRESENT (1-4)
- RCVE UNBLANK (1-4)
- LO PRESENT (5-8)
- RCVE UNBLANK (5-8)

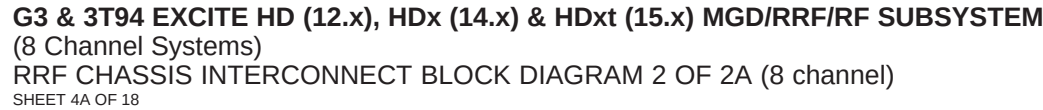


MGD/RRF/RF



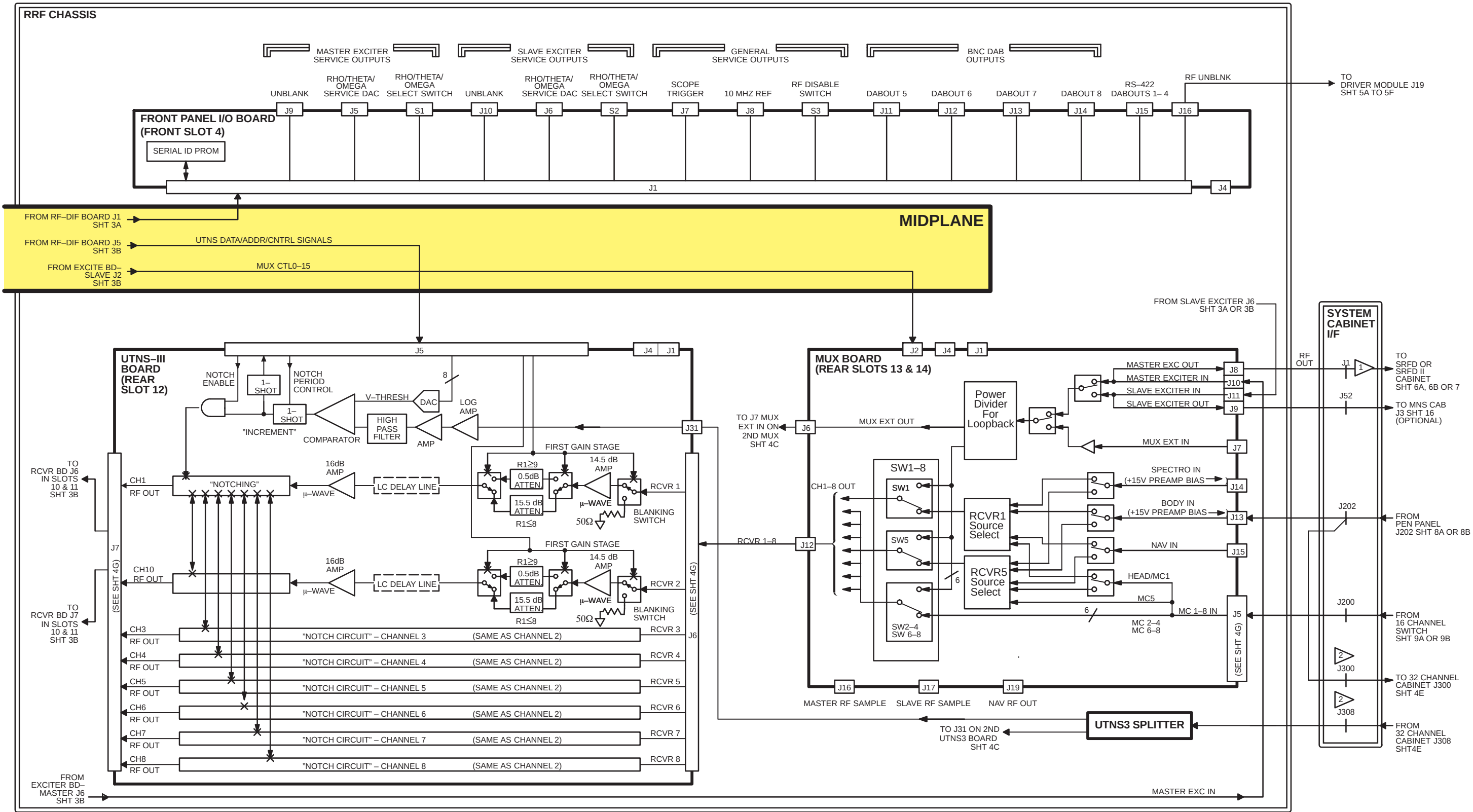


G3 & 3T94 EXCITE HD (12.x), HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
(Channels 1–8 of 32 Channel – Forward Production & Upgrades)
RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 2a OF 4C (32 Channel)
SHEET 4C OF 18


NOTE:
1 J1 not present in RF System Cabinet.
2 32 Channel systems only.

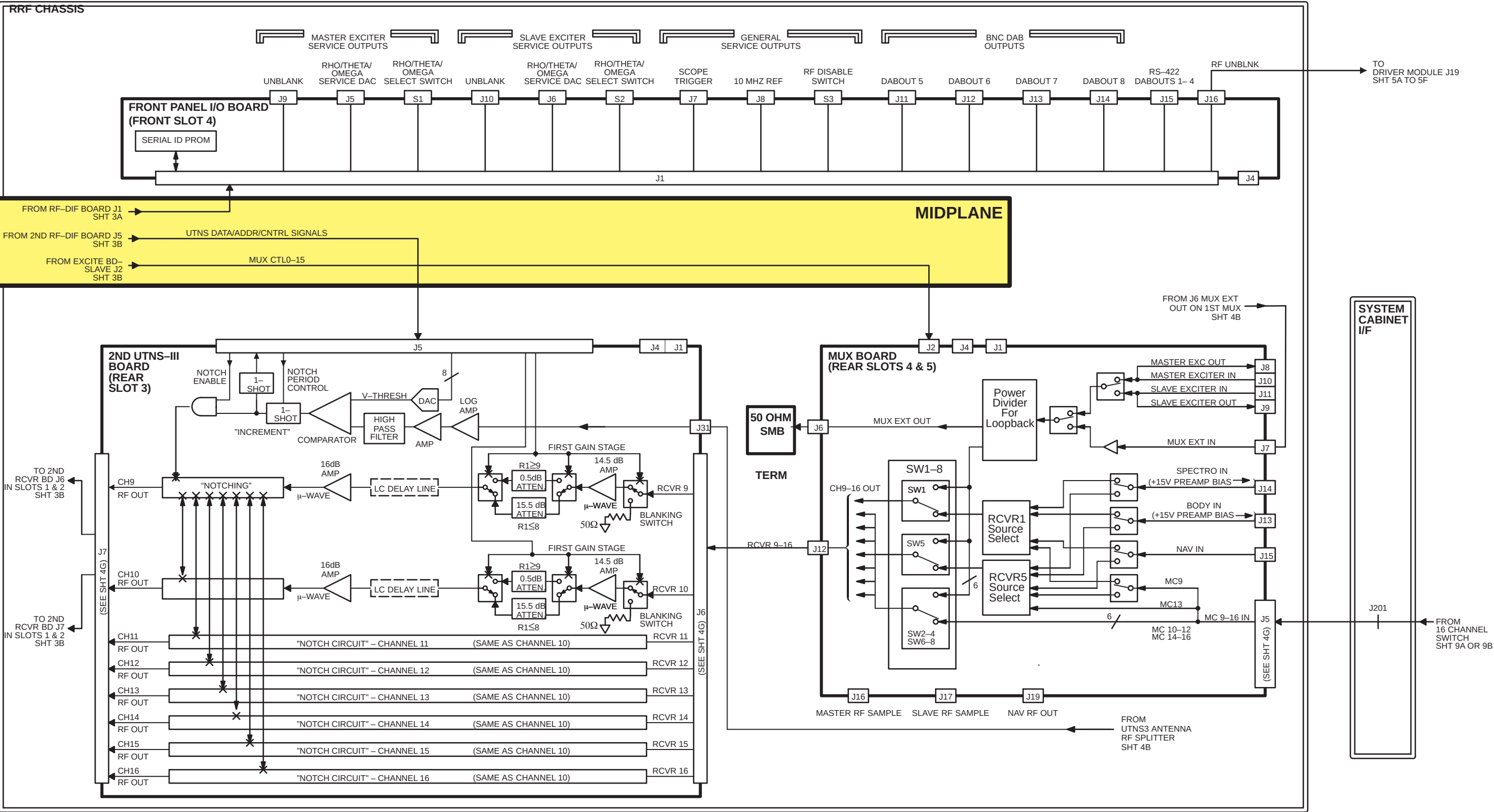


G3 & 3T94 EXCITE HD (12.x) , HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
(Channels 9-16 of 16 or 32 Channel – Forward Production & Upgrades)
RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 2b OF 2B (16 Channel)
RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 2b OF 4C (32 Channel)
SHEET 4D OF 18

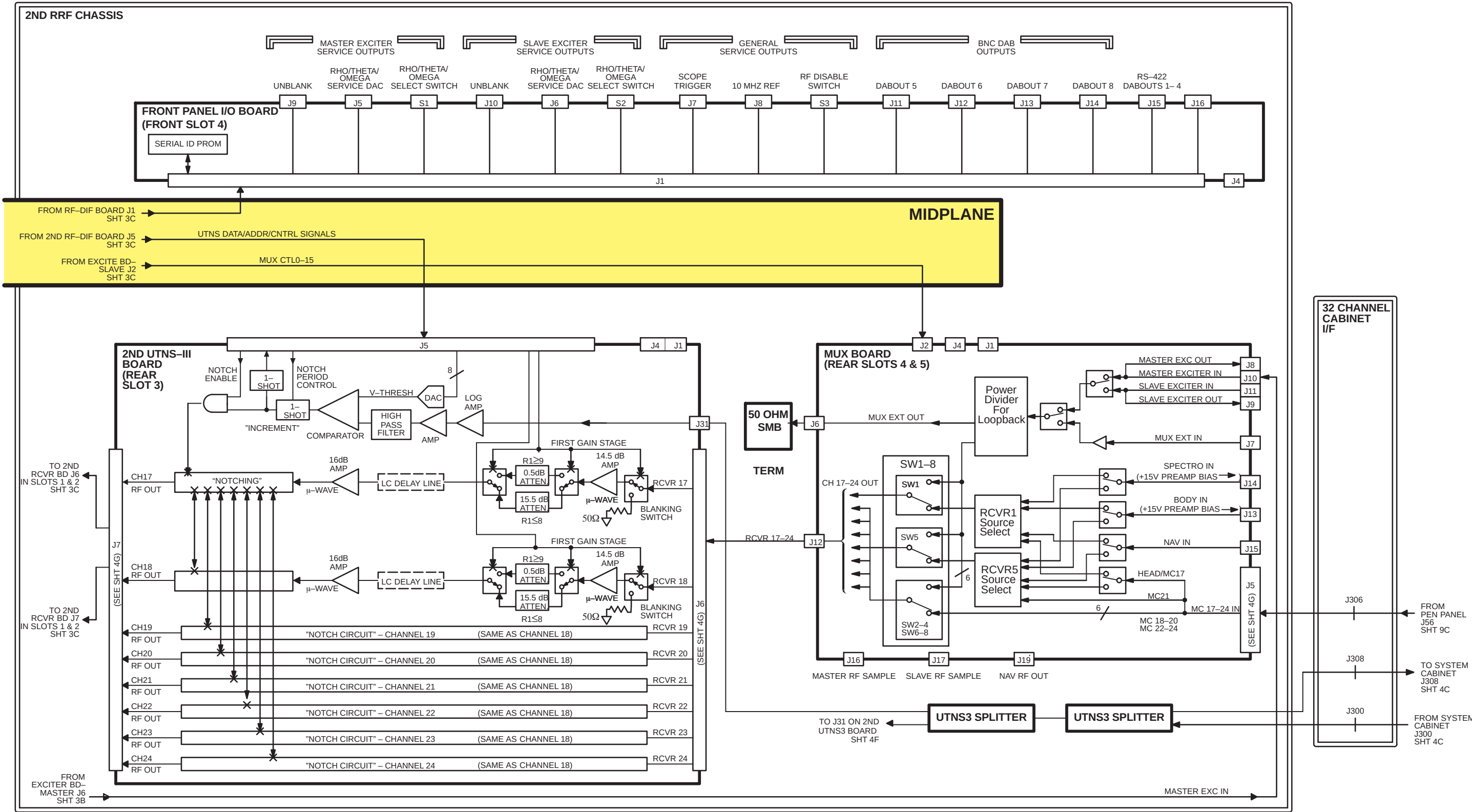

G3 HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
(Channels 17-24 of 32 Channel – Forward Production & Upgrades)
2nd RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 2c OF 4C (32 Channel)
SHEET 4E OF 18

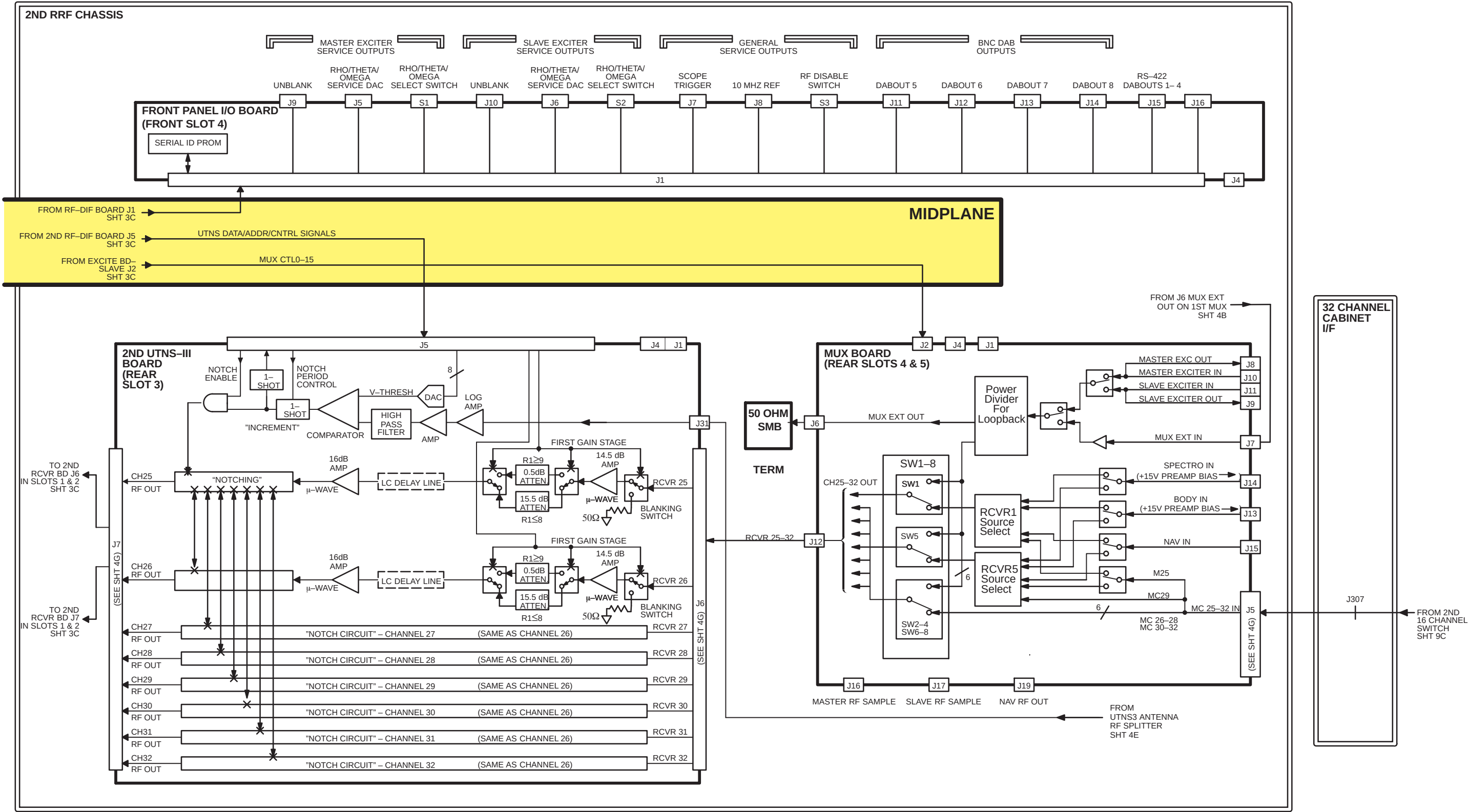

G3 HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
(Channels 25-32 of 32 Channel - Forward Production & Upgrades)
RRF CHASSIS INTERCONNECT BLOCK DIAGRAM 2d OF 4C (32 Channel)
SHEET 4F OF 18

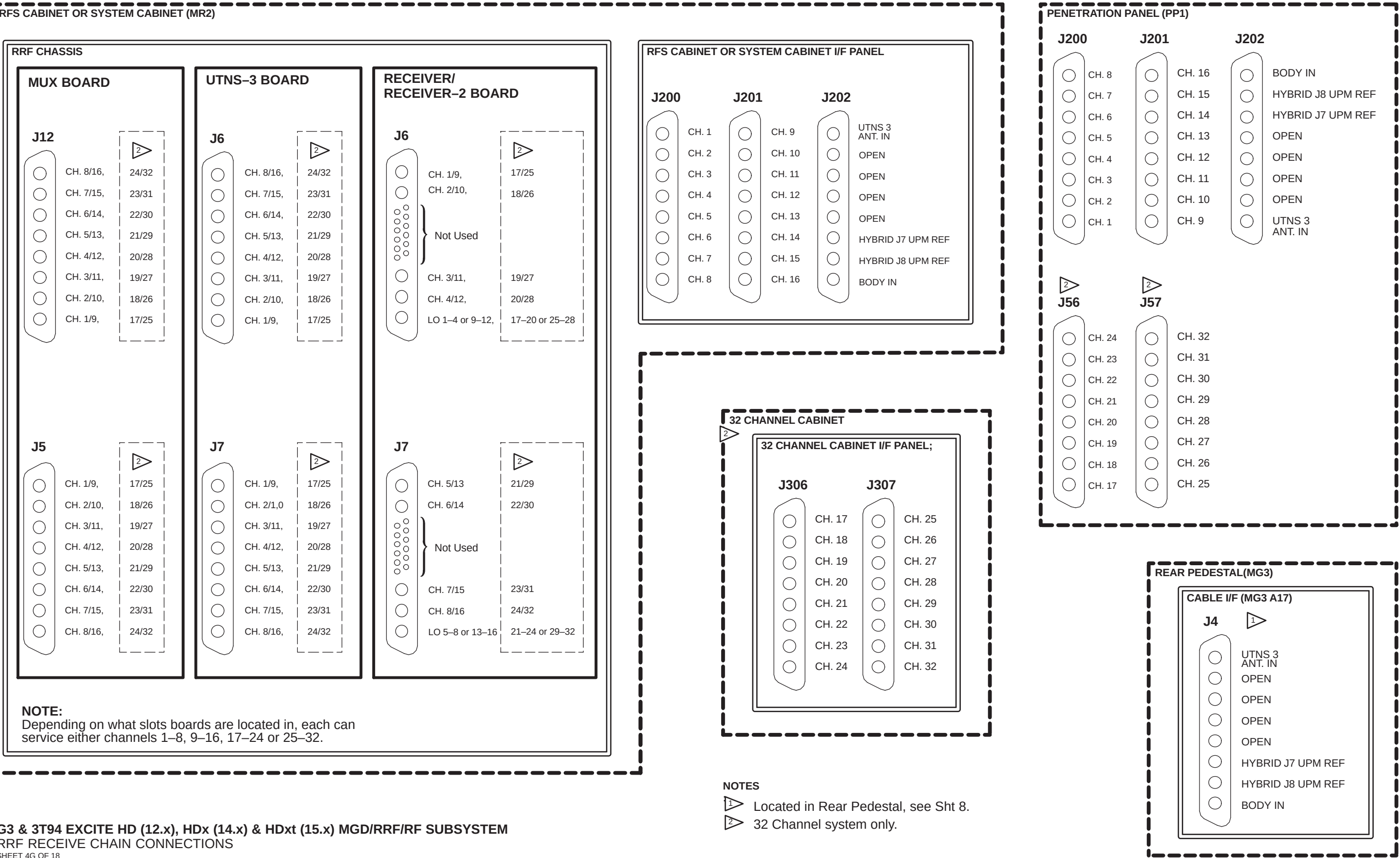

PENETRATION PANEL (PP1)

J200

CH. 8

CH. 7

CH. 6

CH. 5

CH. 4

CH. 3

CH. 2

CH. 1

J201

CH. 16

CH. 15

CH. 14

CH. 13

CH. 12

CH. 11

CH. 10

CH. 9

J202

BODY IN

HYBRID J8 UPM REF

HYBRID J7 UPM REF

OPEN

OPEN

OPEN

OPEN

UTNS 3 ANT. IN

J56

CH. 24

CH. 23

CH. 22

CH. 21

CH. 20

CH. 19

CH. 18

CH. 17

J57

CH. 32

CH. 31

CH. 30

CH. 29

CH. 28

CH. 27

CH. 26

CH. 25

REAR PEDESTAL(MG3)

CABLE I/F (MG3 A17)

J4

UTNS 3 ANT. IN

OPEN

OPEN

OPEN

OPEN

HYBRID J7 UPM REF

HYBRID J8 UPM REF


BODY IN

- NOTES
- 1

 Located in Rear Pedestal, see Sht 8.
- 2


 32 Channel system only.



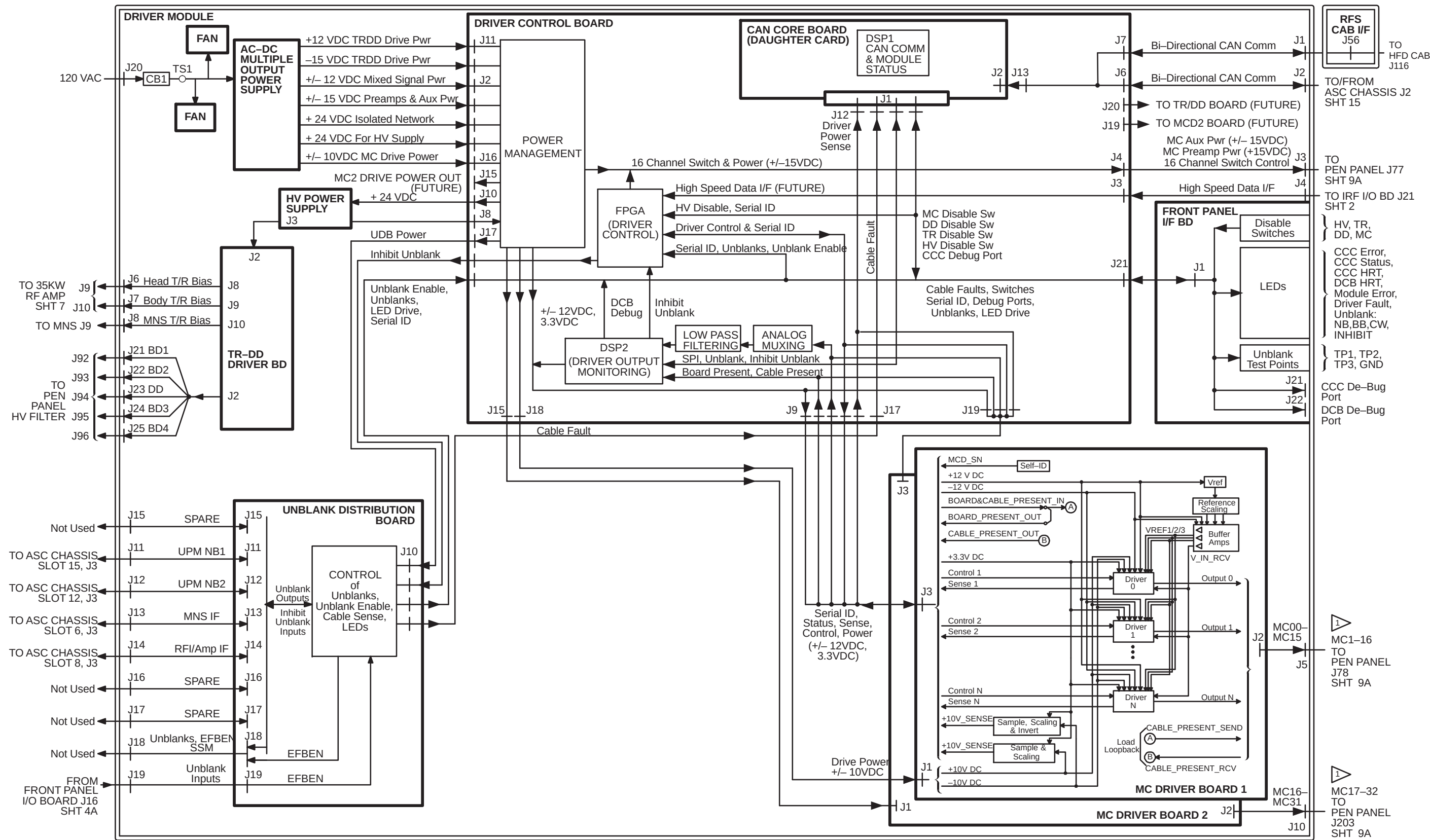


G3 EXCITE HD (12.x) MGD/RRF/RF SUBSYSTEM
DRIVER MODULE (Forward Production, RFS Cabinet, no SSM)
SHEET 5B OF 18

NOTE:
1 MC DRIVER NUMBERING INSIDE THE DRIVER MODULE IS 0-15 AND 16-32. OUTSIDE NUMBERING IS 1-16 AND 17-32.

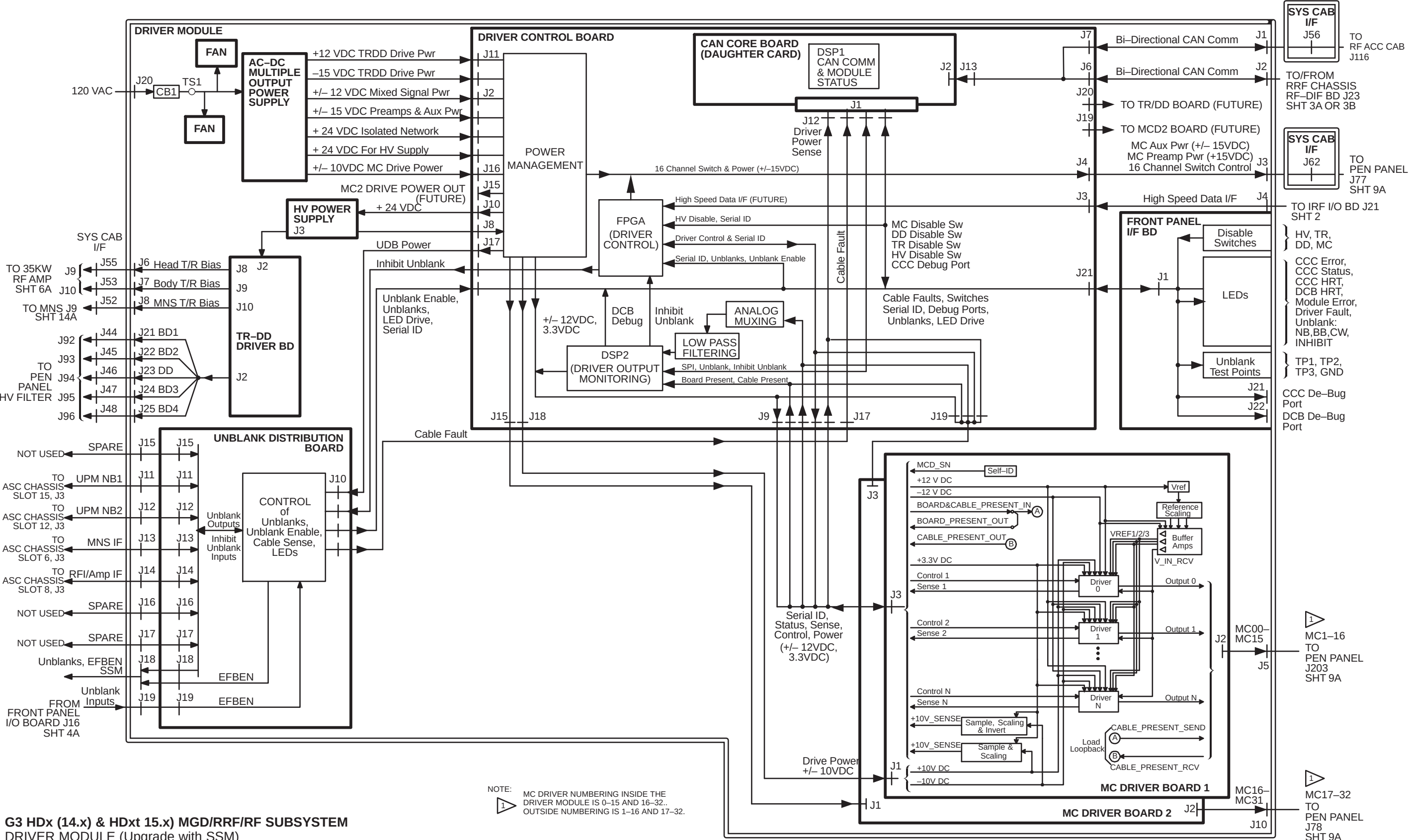


G3 HDx (14.x) & HDxt 15.x) MGD/RRF/RF SUBSYSTEM
DRIVER MODULE (Upgrade with SSM)
SHEET 5C OF 18

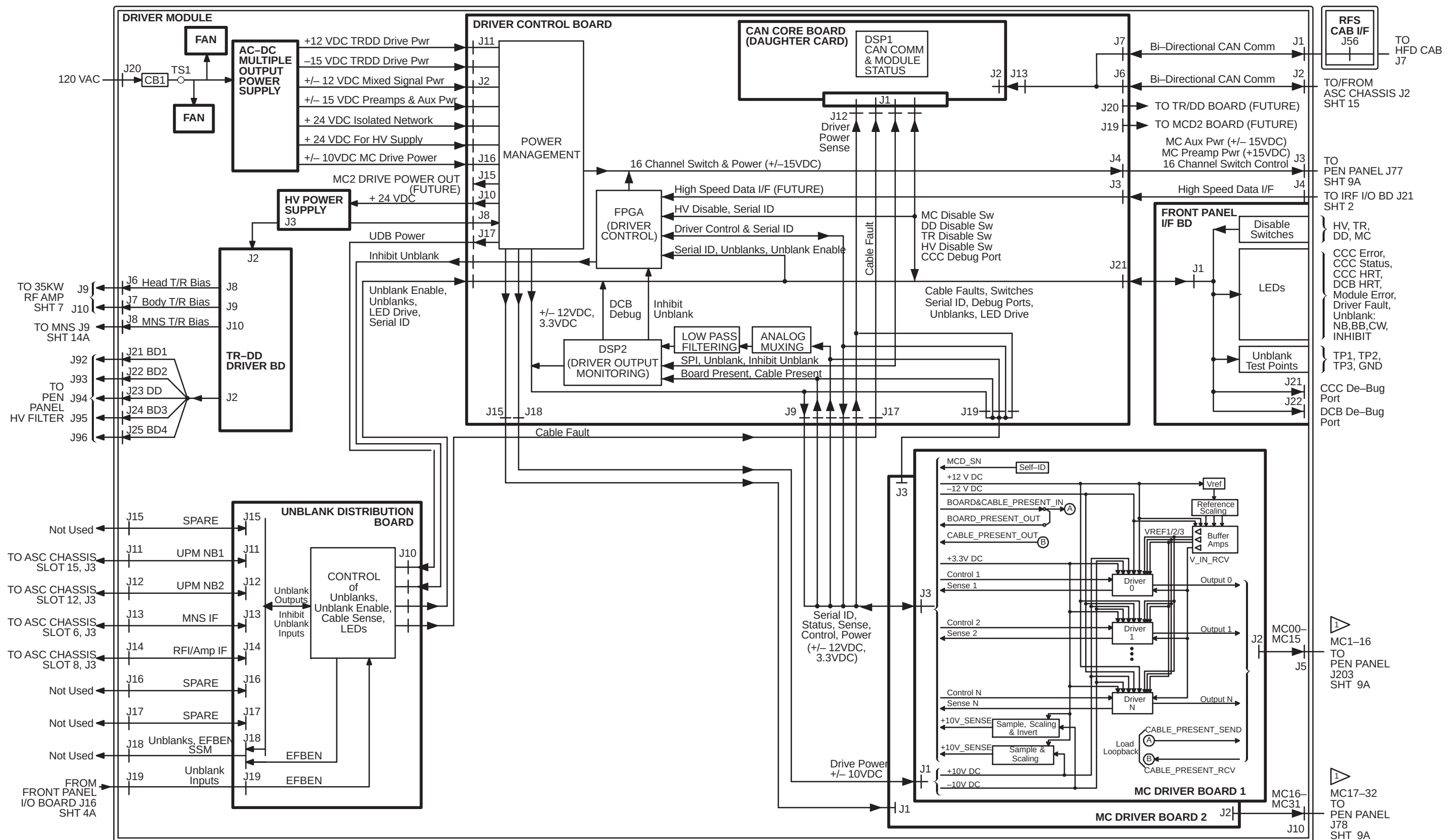


G3 HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
DRIVER MODULE (Forward Production, RFS Cabinet, no SSM)
 SHEET 5D OF 18

NOTE:  MC DRIVER NUMBERING INSIDE THE DRIVER MODULE IS 0-15 AND 16-32.. OUTSIDE NUMBERING IS 1-16 AND 17-32.

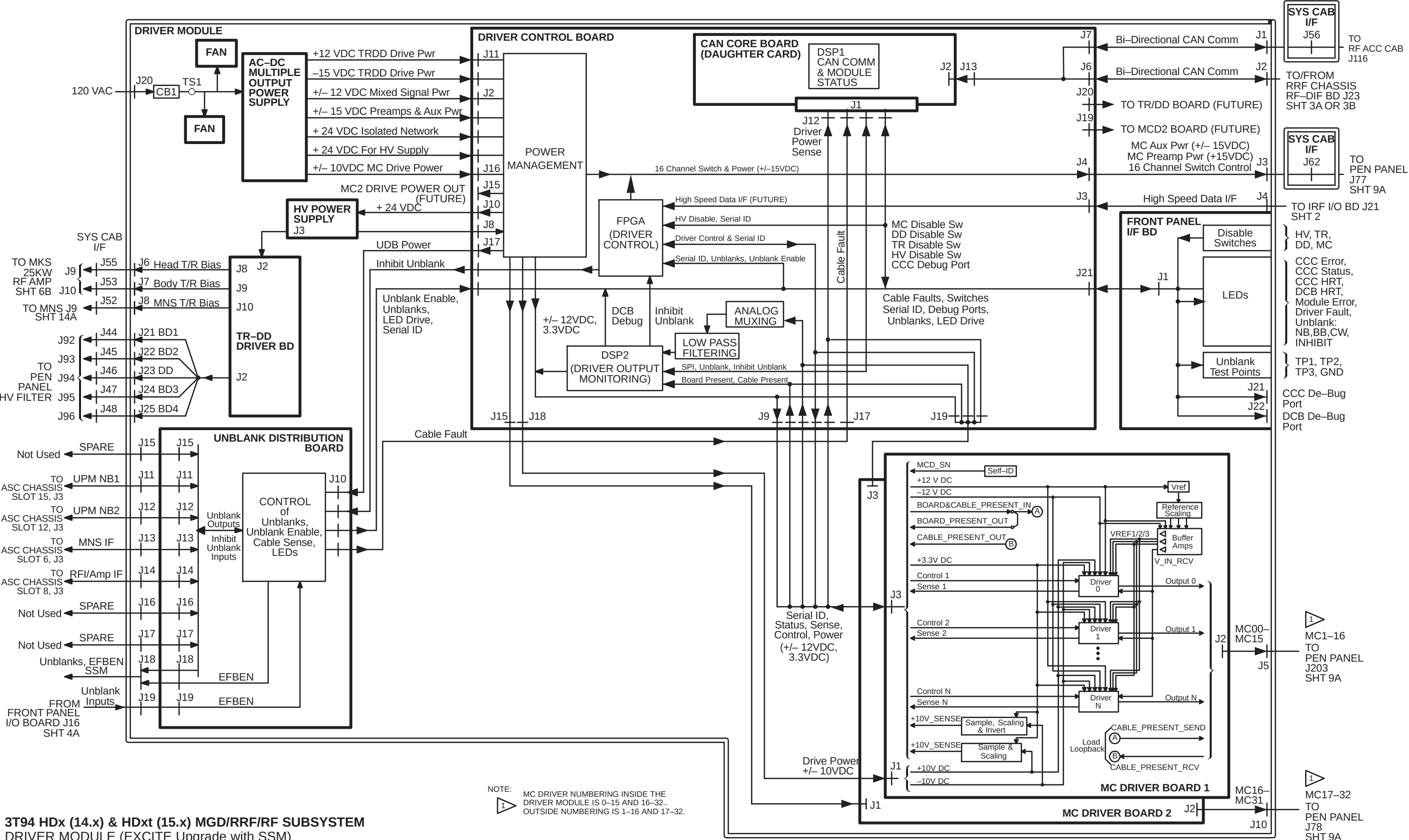


3T94 HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
DRIVER MODULE (EXCITE Upgrade with SSM)
SHEET 5E OF 18

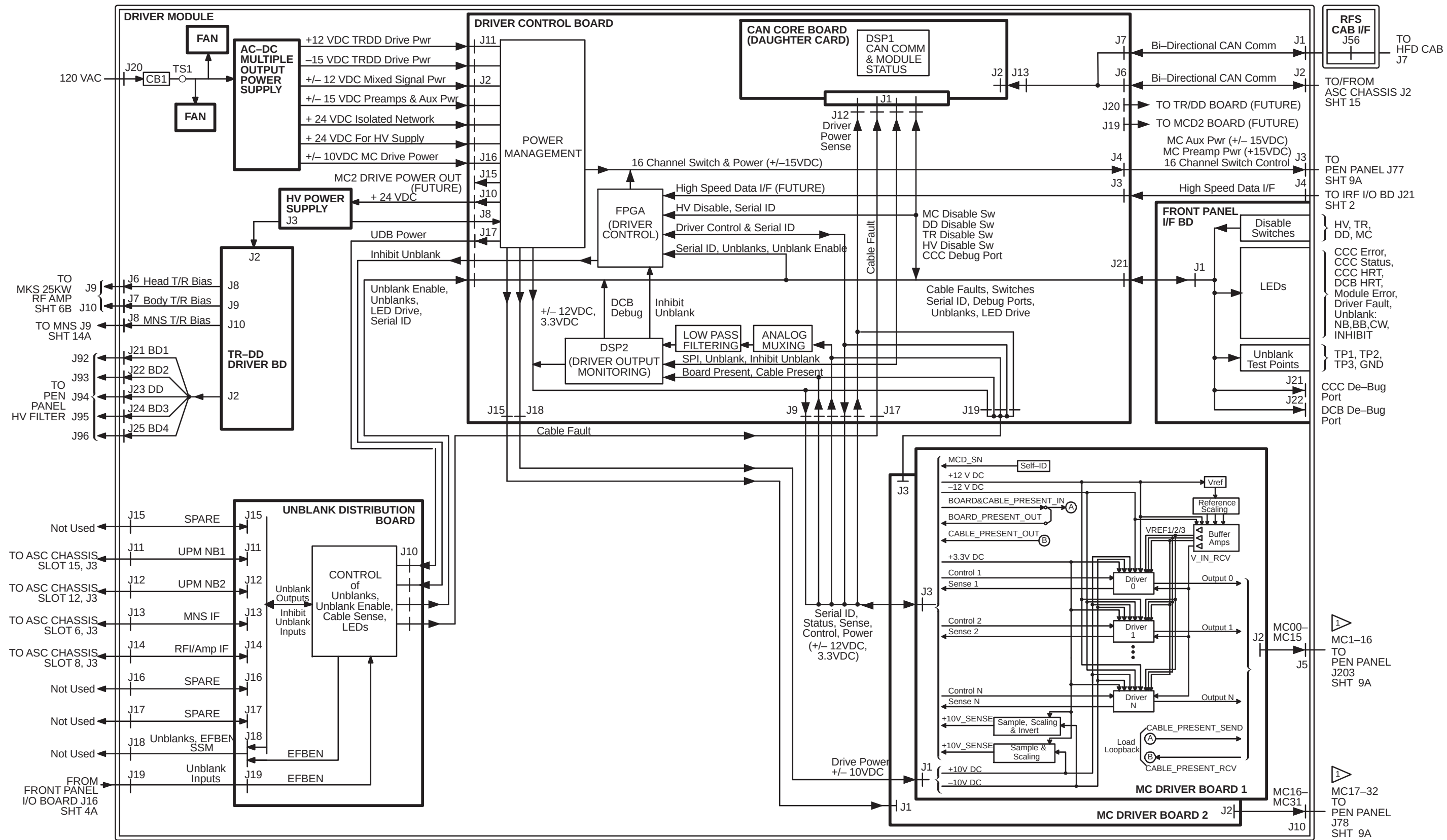


3T94 HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
DRIVER MODULE (Pre-EXCITE UPGRADE, RFS Cabinet, no SSM)
SHEET 5F OF 18


NOTE:
MC DRIVER NUMBERING INSIDE THE
DRIVER MODULE IS 0-15 AND 16-32..
OUTSIDE NUMBERING IS 1-16 AND 17-32.




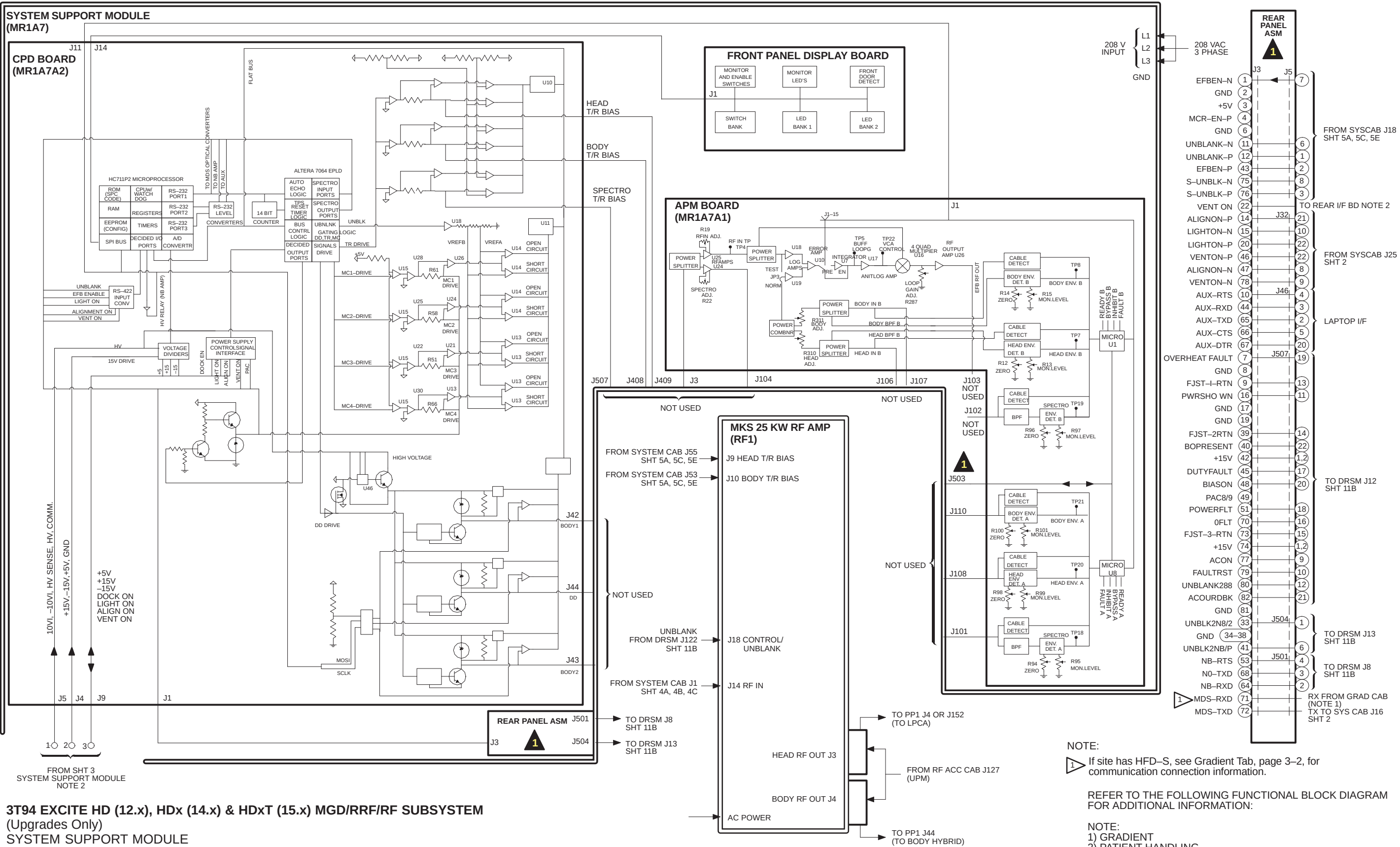

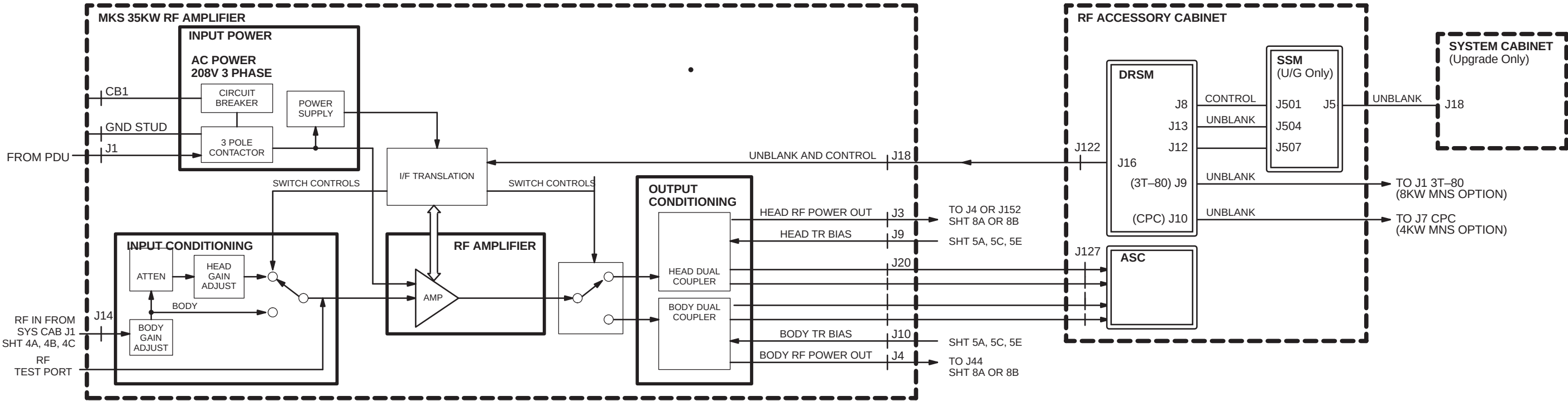



OR

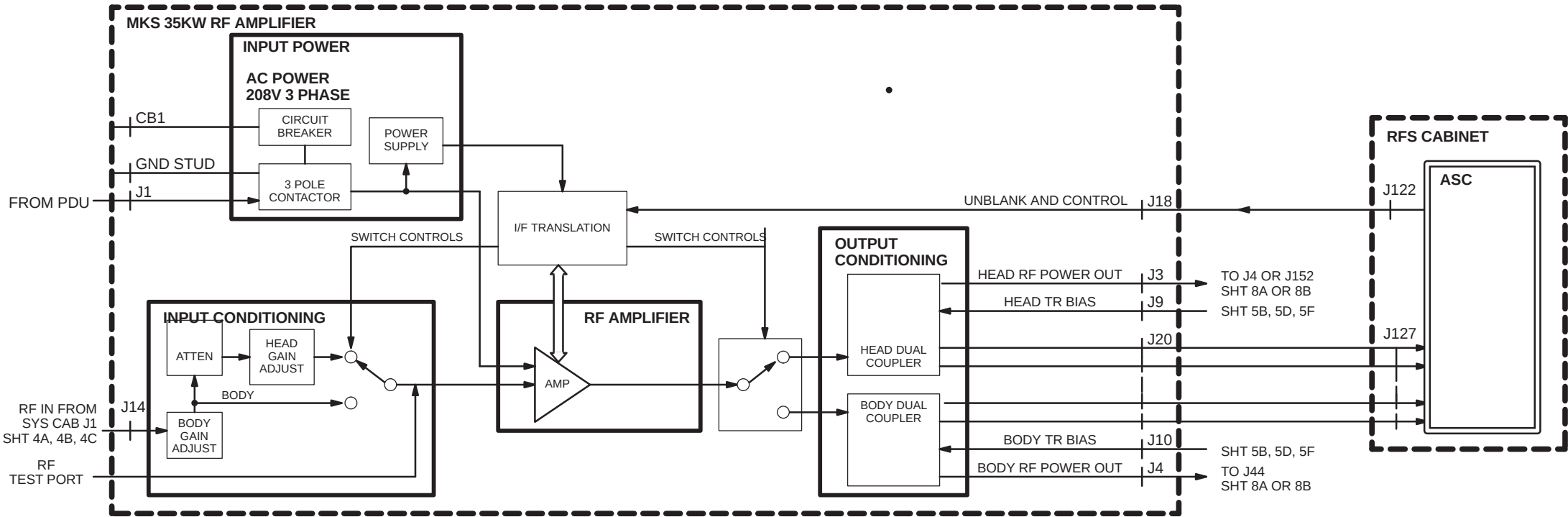

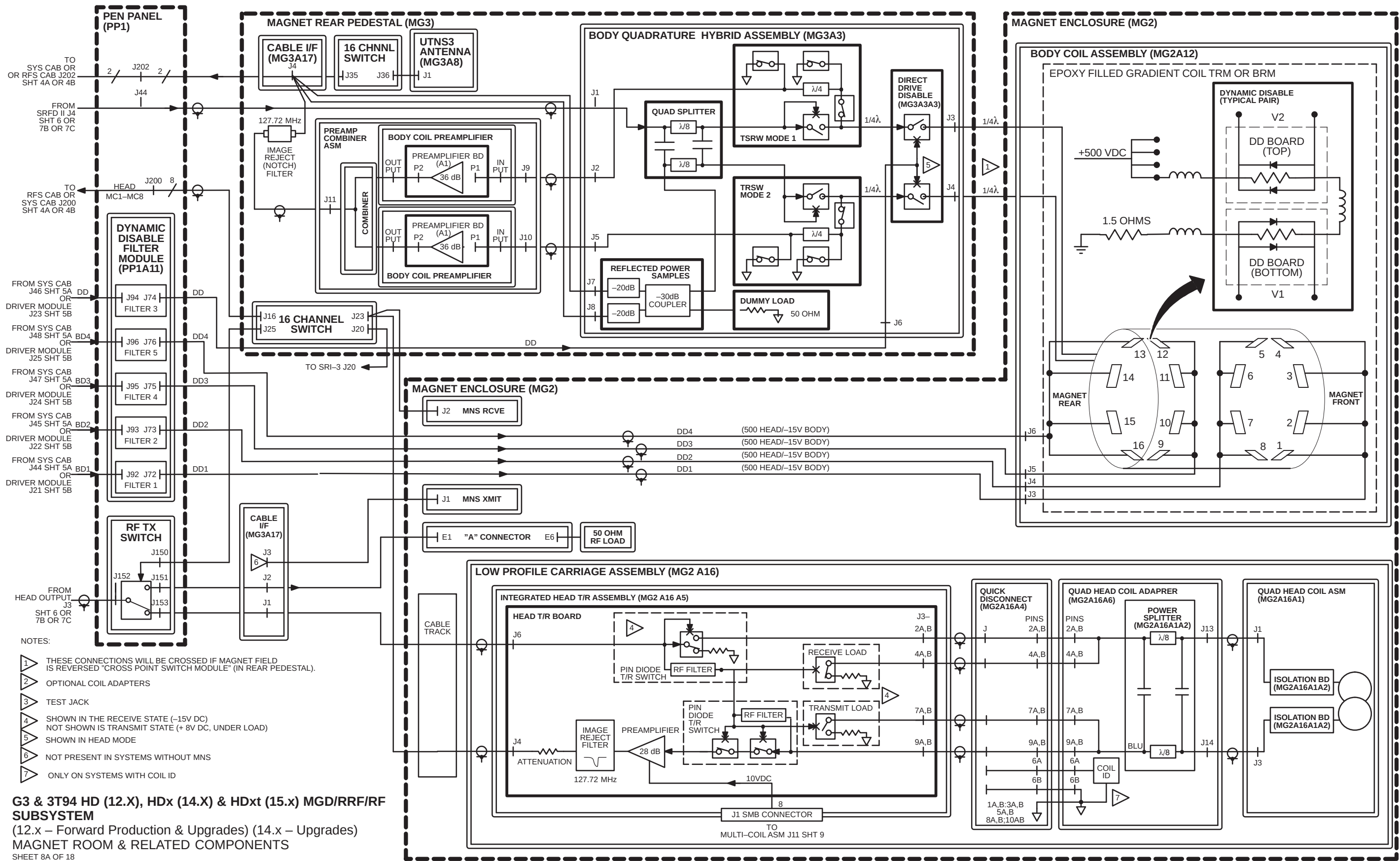


G3 & 3T94 HD (12.X), HDx (14.X) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
(12.x – Forward Production & Upgrades) (14.x – Upgrades)
MAGNET ROOM & RELATED COMPONENTS
SHEET 8A OF 18

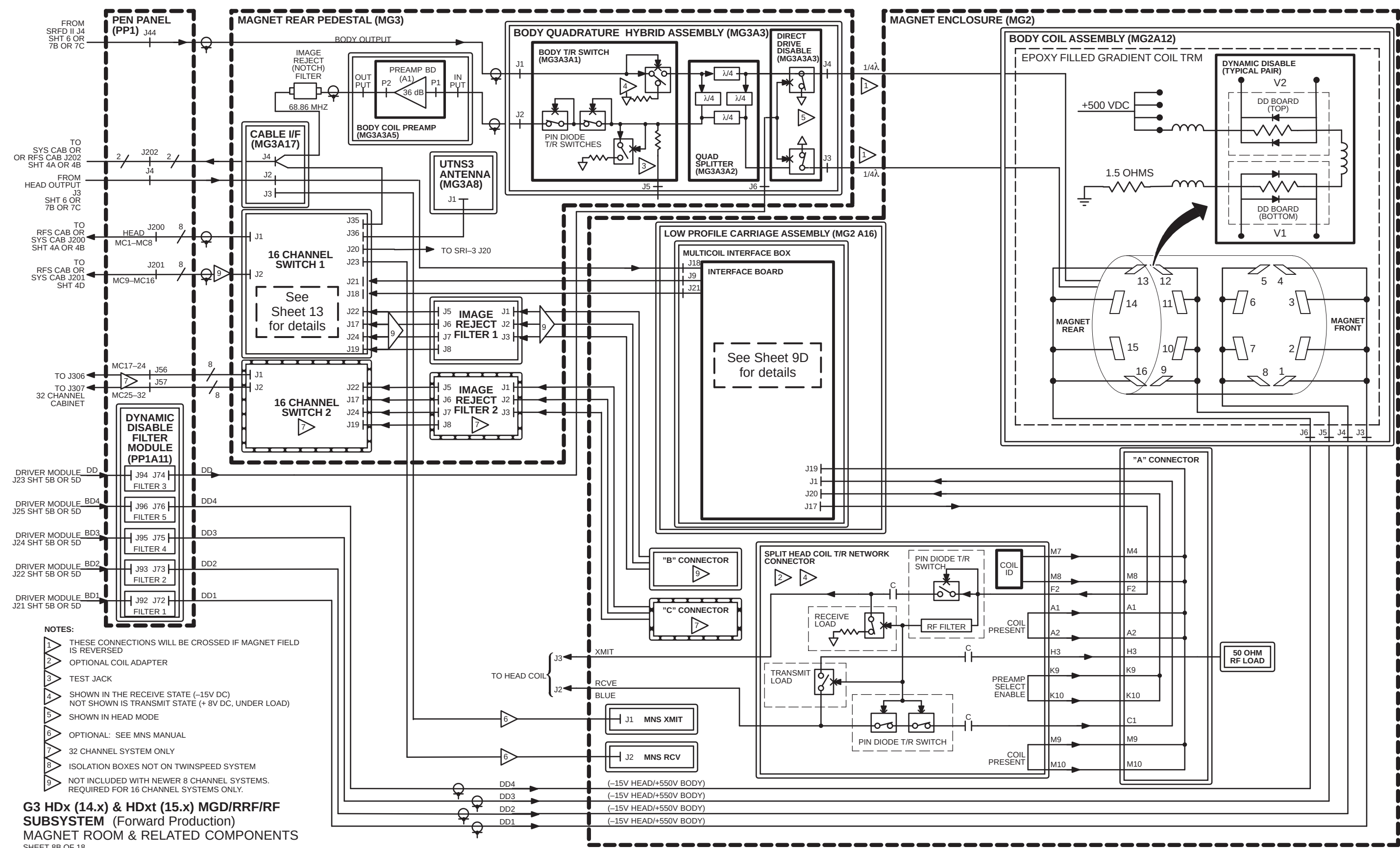

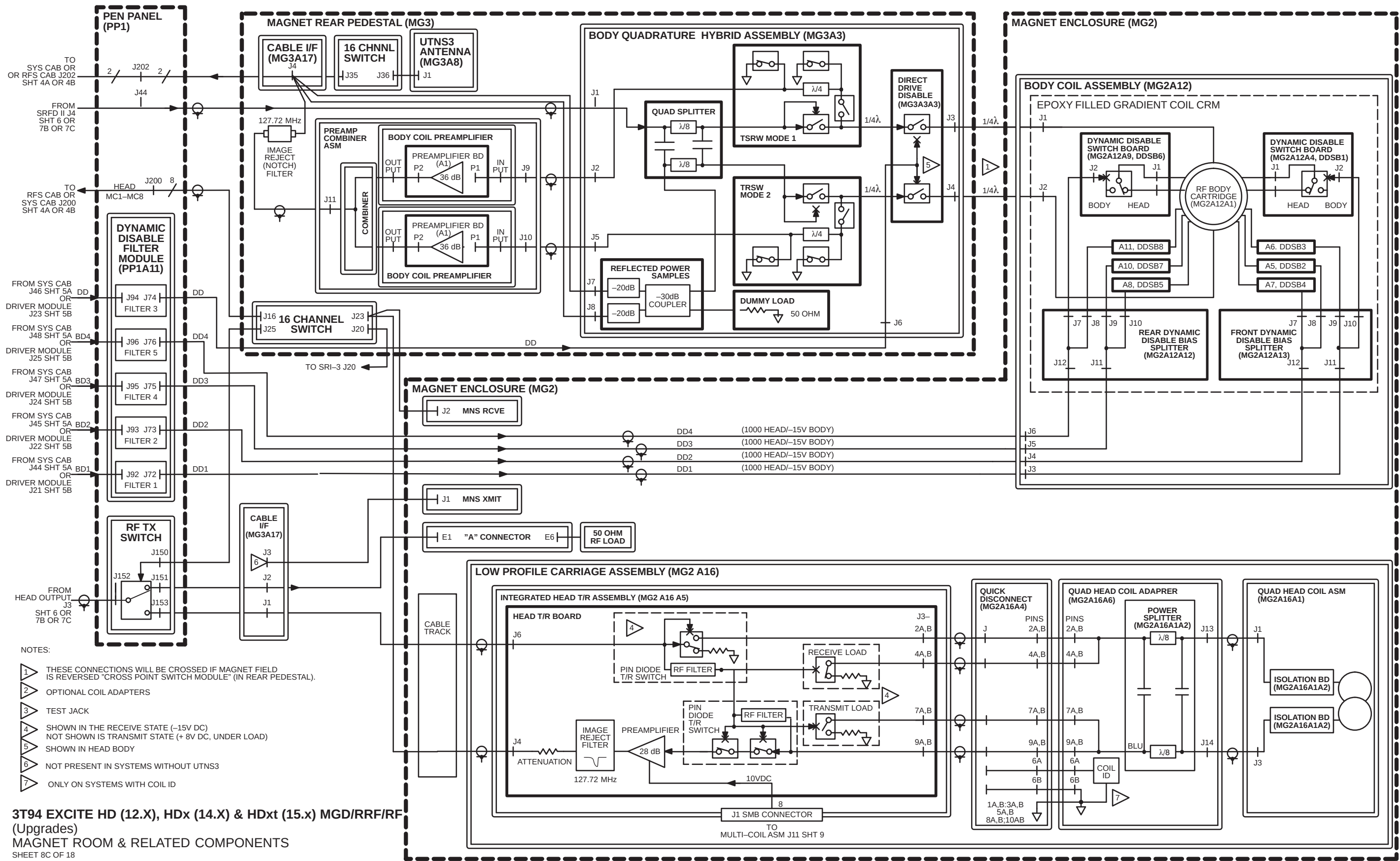





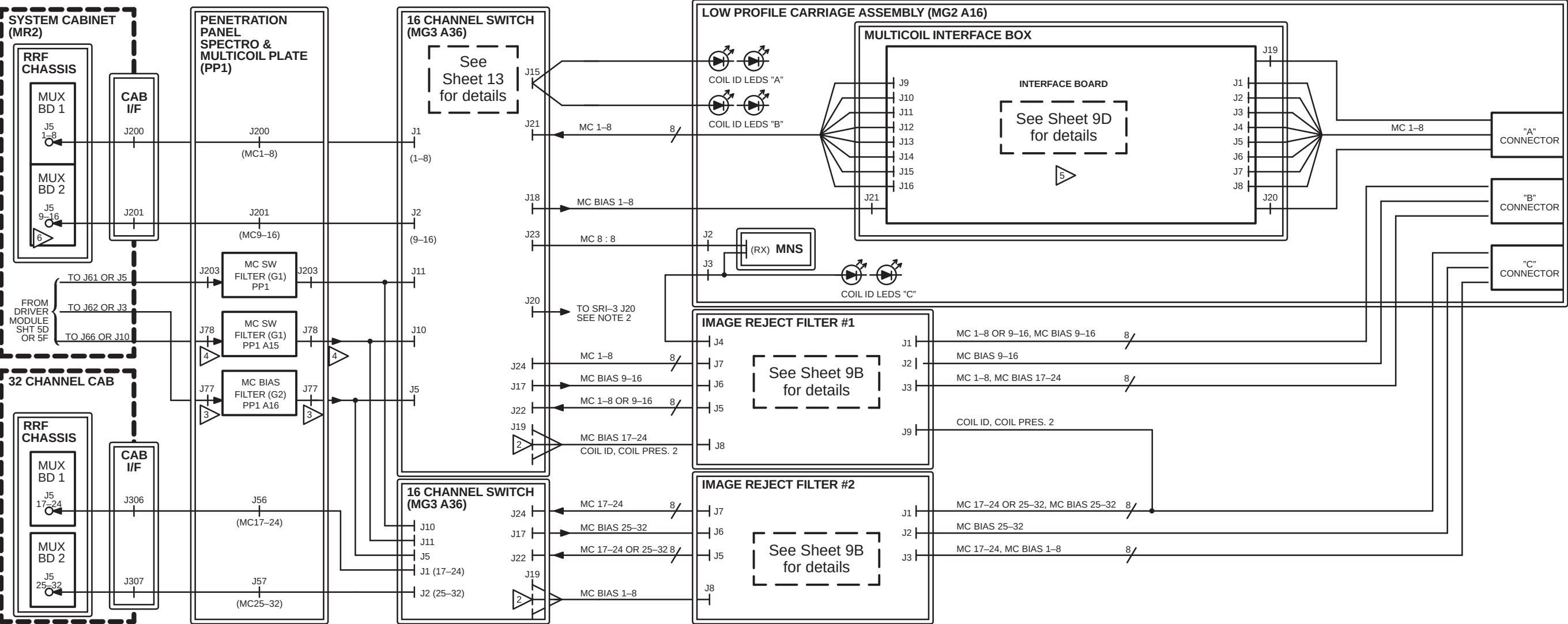

REFER TO THE FOLLOWING FUNCTIONAL BLOCK DIAGRAM FOR ADDITIONAL INFORMATION:

NOTE:

1) OPERATOR WORKSPACE
2) PATIENT HANDLING

NOTES:

1

G3 HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
32 CHANNEL MULTI-COIL

SHEET 9C OF 18

2

J8 & J19

PIN(S) SIGNALS

1,14 -15V
2,15 +15V
3, 4, 23-25 . . . GND
5 MCD8
6 MCD4
7 MCD7
8 MCD3
9 MCD6
10 MCD2
11 MCD5
12 MCD1
13 CTL_PREAMP (NOTE 1)
16 RED_LED
17 YEL_LED
18 GRN_LED
19 COIL PRESENT #1
20 COIL PRESENT #2
21 COIL ID #1
22 COIL ID #2

NOTE 1: PIN 13 DEFAULT LOW.
LOW: PREAMP POWER ON.
HIGH: PREAMP POWER OFF.

3

J62 OR J3 & J77

PIN SIGNAL NAME

1 XPT_DATA_POS
2 XPT_STROBE_POS
3 XPT_OE_N_POS
4 AGND
5 AGND
6 +15V (NOTE 1)
7 +15V (NOTE 1)
8 +15V (NOTE 1)
9 XPT_DATA_NEG
10 XPT_STROBE_NEG
11 XPT_OE_N_NEG
12 AGND
13 AGND
14 -15V (MC COILS)
15 -15V (MC COILS)

NOTE 1: ANALOG POWER FOR
PREAMPS (MC & HEAD),
CROSSPOINT & MC COILS.

4

J10 OR J66 & J78

PIN(S) SIGNAL

1-8 MCD1-8 (NOTE 1)
9 MC CABLE PRESENT_O IN (NOTE 2)
10,11 NO CONNECT
12-19 MCD9-16 (NOTE 3)
20-27 GND
29 MC CABLE PRESENT_1 IN (NOTE 2)
30-37 GND

NOTE 1: MULTICOIL BIAS (+3V,+5V,+7V,-5V).
NOTE 2: THESE TWO PINS CONNECTED TOGETHER
ON J10 FOR DIAGNOSTIC LOOPBACK.
NOTE 3: NOT USED/CONNECTED AT J10 INSIDE
CROSSPOINT SWITCH ASSEMBLY.


5


TESTPOINT . FUNCTION

TP1 SYSTEM PREAMP ENABLE (+15V OR 0V)
TP2 PREAMP SELECT 1
TP3 PREAMP SELECT 2
TP4 8.5V_A
TP5 +10V_B CNTRL
TP6 +15V
TP7 +5.1V
TP8 +15V
TP9 +15V
TP10 8.5V_B
TP11 8.5V_C
TP12 -15V
TP13 -8.5V

6

Only 16 Channel Systems have a second MUX board. J201 not used in 8 and 4 Channel systems.

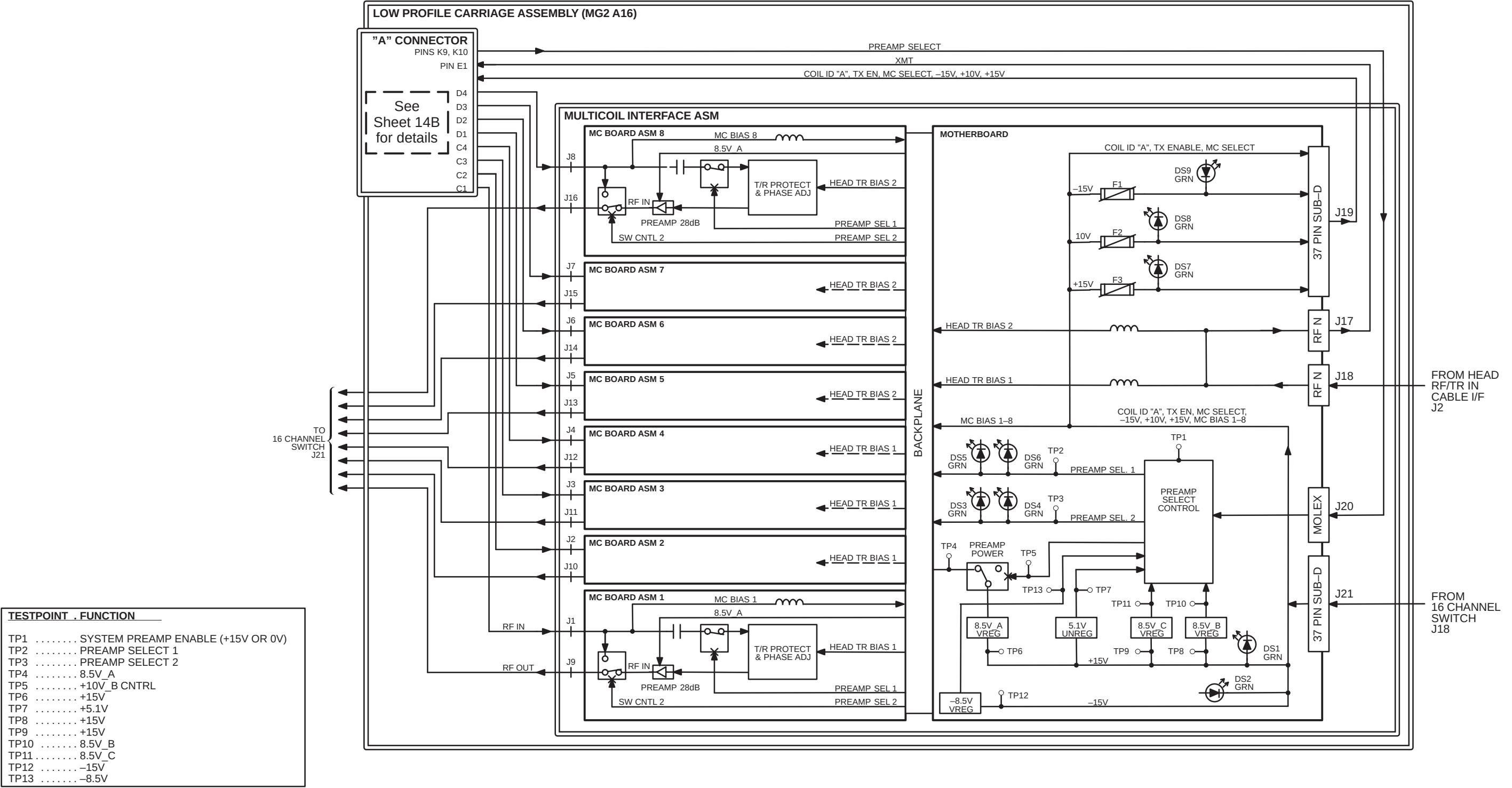

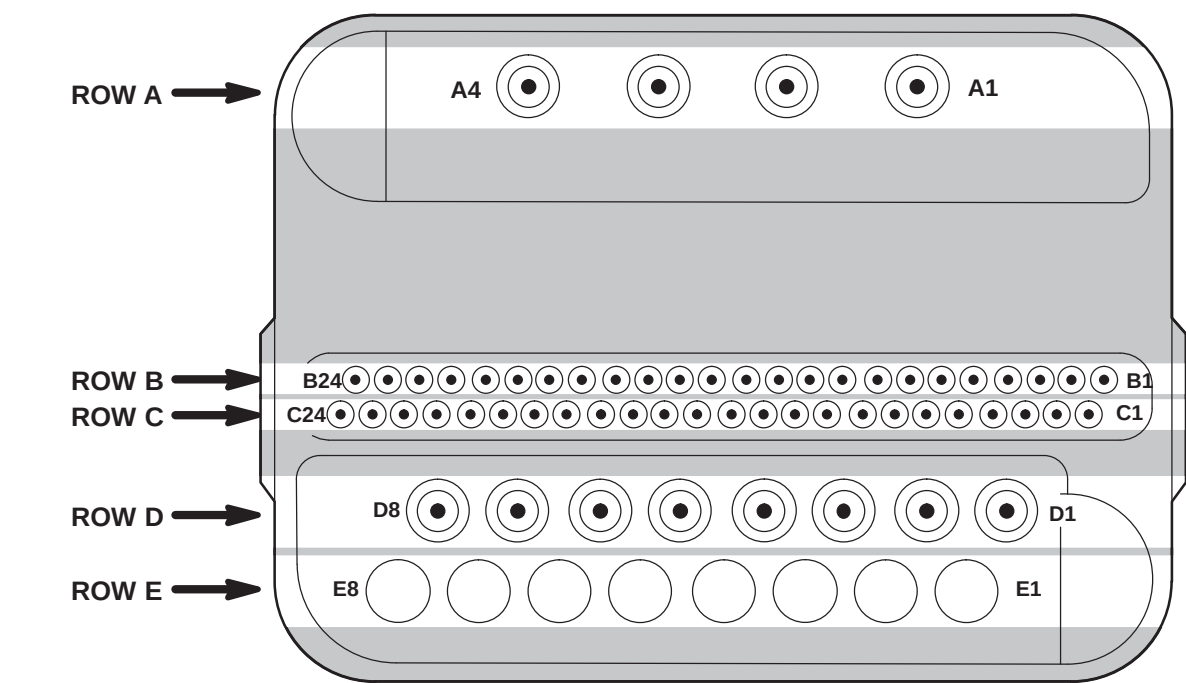

As viewed from the front of the Hypertronic connector

HYPERTRONICS CONNECTOR
(ROW A PINOUTS)
HEAD AND SURFACE COIL TO SYSTEM INTERCONNECT

Pin Out		Description	Style	PreAmp on Coil Connection Requirements	System PreAmp Connection Requirements
A	1	Head Coil Receive w/Transmit Load	Co–Axial	N/A	
A	2	N/A		N/A	
A	3	N/A		N/A	
A	4	Head Coil Transmit w/Transmit Load	Co–Axial	N/A	

G3 & 3T94 EXCITE HD (12.x), HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
HYPERTRONIC CONNECTOR PINOUTS
(12.x – Forward Production & Upgrades)
(14.x – Upgrades)
SHEET 9E OF 18


NOTE: For Pin–outs B10 thru B17 below, MC Bias is 17–24 on 14.x systems.


HYPERTRONICS CONNECTOR
(ROW B PINOUTS)
DC, COIL ID, AND MULTICOIL BIAS TO SYSTEM INTERCONNECT

Pin Out		Description	Style	PreAmp on Coil Connection Requirements	System PreAmp Connection Requirements
B	1	+ 15 V DC	Pin		
B	2	Ground	Pin	Ground	Ground
B	3	+10 V_BPreAmp DC	Pin	+10 Volt 1A for PreAmp coil.	N/C
B	4	Ground	Pin	Ground	Ground
B	5	Spare	Pin		
B	6	Spare	Pin		
B	7	Spare	Pin		
B	8	Coil ID	Pin		
B	9	Ground	Pin		
B	10	T/R Switch Bias 1 (17)	Pin	T/R Switch Bias	T/R Switch Bias
B	11	T/R Switch Bias 2 (18)	Pin	T/R Switch Bias	T/R Switch Bias
B	12	T/R Switch Bias 3 (19)	Pin	T/R Switch Bias	T/R Switch Bias
B	13	T/R Switch Bias 4 (20)	Pin	T/R Switch Bias	T/R Switch Bias
B	14	T/R Switch Bias 5 (21)	Pin	T/R Switch Bias	T/R Switch Bias
B	15	T/R Switch Bias 6 (22)	Pin	T/R Switch Bias	T/R Switch Bias
B	16	T/R Switch Bias 7 (23)	Pin	T/R Switch Bias	T/R Switch Bias
B	17	T/R Switch Bias 8 (24)	Pin	T/R Switch Bias	T/R Switch Bias
B	18	Integrated Preamp (System) Switch Control	Pin	For Preamp on Coil: Connect B18 with C17	For Non–Preamp on Coil: Leave this pin as Non–Connect
B	19	Spare	Pin		
B	20	Coil Present	Pin	Must connect B20 to C6	Must connect B20 to C6
B	21	Spare	Pin		
B	22	+15 VDC	Pin	1A max	1A max
B	23	Spare	Pin		
B	24	Spare	Pin		


GE Healthcare																																																																																																																																																											
REV 12																																																																																																																																																											
Signal HD, HDx & HDxt 3.0T (12.x, 14.x, 15.x & 16.x) Block Diagrams & Supplemental Schematics																																																																																																																																																											
DIRECTION 5133320																																																																																																																																																											
HYPERTRONICS CONNECTOR (ROW C PINOUTS) DC, COIL ID, AND MULTICOIL TO SYSTEM INTERCONNECT																																																																																																																																																											
HYPERTRONICS CONNECTOR (ROW C PINOUTS) DC, COIL ID, AND MULTICOIL TO SYSTEM INTERCONNECT																																																																																																																																																											
<table><tr><th colspan="2">Pin Out</th><th>Description</th><th>Style</th><th>PreAmp on Coil Connection Requirements</th><th>System PreAmp Connection Requirements</th></tr><tr><td>C</td><td>1</td><td>Spare</td><td>Pin</td><td></td><td></td></tr><tr><td>C</td><td>2</td><td>Spare</td><td>Pin</td><td></td><td></td></tr><tr><td>C</td><td>3</td><td>Spare</td><td>Pin</td><td></td><td></td></tr><tr><td>C</td><td>4</td><td>Spare</td><td>Pin</td><td></td><td></td></tr><tr><td>C</td><td>5</td><td>Spare</td><td>Pin</td><td></td><td></td></tr><tr><td>C</td><td>6</td><td>Coil Present Return</td><td>Pin</td><td>Must connect B20 to C6</td><td>Must connect B20 to C6</td></tr><tr><td>C</td><td>7</td><td>Coil ID Return</td><td>Pin</td><td></td><td></td></tr><tr><td>C</td><td>8</td><td>Ground</td><td>Pin</td><td>Ground</td><td>Ground</td></tr><tr><td>C</td><td>9</td><td>Reserved</td><td>Pin</td><td>Reserved</td><td>Reserved</td></tr><tr><td>C</td><td>10</td><td>Reserved</td><td>Pin</td><td>Reserved</td><td>Reserved</td></tr><tr><td>C</td><td>11</td><td>Reserved</td><td>Pin</td><td>Reserved</td><td>Reserved</td></tr><tr><td>C</td><td>12</td><td>Reserved</td><td>Pin</td><td>Reserved</td><td>Reserved</td></tr><tr><td>C</td><td>13</td><td>Reserved</td><td>Pin</td><td>Reserved</td><td>Reserved</td></tr><tr><td>C</td><td>14</td><td>Reserved</td><td>Pin</td><td>Reserved</td><td>Reserved</td></tr><tr><td>C</td><td>15</td><td>Reserved</td><td>Pin</td><td>Reserved</td><td>Reserved</td></tr><tr><td>C</td><td>16</td><td>Reserved</td><td>Pin</td><td>Reserved</td><td>Reserved</td></tr><tr><td>C</td><td>17</td><td>Inegrated Preamp Switch Control RTN</td><td>Pin</td><td>For Preamp on Coil: Connect B18 with C17</td><td>For Non–Preamp on coil, leave this pin as Non–Connect</td></tr><tr><td>C</td><td>18</td><td>Ground</td><td>Pin</td><td>Ground</td><td>Ground</td></tr><tr><td>C</td><td>19</td><td>Spare</td><td>Pin</td><td></td><td></td></tr><tr><td>C</td><td>20</td><td>Ground</td><td>Pin</td><td>Ground</td><td>Ground</td></tr><tr><td>C</td><td>21</td><td>+10V_BPreAmp DC</td><td>Pin</td><td>+10 Volt 1A for Preamp Coil</td><td>N/C</td></tr><tr><td>C</td><td>22</td><td>Ground</td><td>Pin</td><td>Ground</td><td>Ground</td></tr><tr><td>C</td><td>23</td><td>Spare</td><td>Pin</td><td></td><td></td></tr><tr><td>C</td><td>24</td><td>Spare</td><td>Pin</td><td></td><td></td></tr></table>						Pin Out		Description	Style	PreAmp on Coil Connection Requirements	System PreAmp Connection Requirements	C	1	Spare	Pin			C	2	Spare	Pin			C	3	Spare	Pin			C	4	Spare	Pin			C	5	Spare	Pin			C	6	Coil Present Return	Pin	Must connect B20 to C6	Must connect B20 to C6	C	7	Coil ID Return	Pin			C	8	Ground	Pin	Ground	Ground	C	9	Reserved	Pin	Reserved	Reserved	C	10	Reserved	Pin	Reserved	Reserved	C	11	Reserved	Pin	Reserved	Reserved	C	12	Reserved	Pin	Reserved	Reserved	C	13	Reserved	Pin	Reserved	Reserved	C	14	Reserved	Pin	Reserved	Reserved	C	15	Reserved	Pin	Reserved	Reserved	C	16	Reserved	Pin	Reserved	Reserved	C	17	Inegrated Preamp Switch Control RTN	Pin	For Preamp on Coil: Connect B18 with C17	For Non–Preamp on coil, leave this pin as Non–Connect	C	18	Ground	Pin	Ground	Ground	C	19	Spare	Pin			C	20	Ground	Pin	Ground	Ground	C	21	+10V_BPreAmp DC	Pin	+10 Volt 1A for Preamp Coil	N/C	C	22	Ground	Pin	Ground	Ground	C	23	Spare	Pin			C	24	Spare	Pin		
Pin Out		Description	Style	PreAmp on Coil Connection Requirements	System PreAmp Connection Requirements																																																																																																																																																						
C	1	Spare	Pin																																																																																																																																																								
C	2	Spare	Pin																																																																																																																																																								
C	3	Spare	Pin																																																																																																																																																								
C	4	Spare	Pin																																																																																																																																																								
C	5	Spare	Pin																																																																																																																																																								
C	6	Coil Present Return	Pin	Must connect B20 to C6	Must connect B20 to C6																																																																																																																																																						
C	7	Coil ID Return	Pin																																																																																																																																																								
C	8	Ground	Pin	Ground	Ground																																																																																																																																																						
C	9	Reserved	Pin	Reserved	Reserved																																																																																																																																																						
C	10	Reserved	Pin	Reserved	Reserved																																																																																																																																																						
C	11	Reserved	Pin	Reserved	Reserved																																																																																																																																																						
C	12	Reserved	Pin	Reserved	Reserved																																																																																																																																																						
C	13	Reserved	Pin	Reserved	Reserved																																																																																																																																																						
C	14	Reserved	Pin	Reserved	Reserved																																																																																																																																																						
C	15	Reserved	Pin	Reserved	Reserved																																																																																																																																																						
C	16	Reserved	Pin	Reserved	Reserved																																																																																																																																																						
C	17	Inegrated Preamp Switch Control RTN	Pin	For Preamp on Coil: Connect B18 with C17	For Non–Preamp on coil, leave this pin as Non–Connect																																																																																																																																																						
C	18	Ground	Pin	Ground	Ground																																																																																																																																																						
C	19	Spare	Pin																																																																																																																																																								
C	20	Ground	Pin	Ground	Ground																																																																																																																																																						
C	21	+10V_BPreAmp DC	Pin	+10 Volt 1A for Preamp Coil	N/C																																																																																																																																																						
C	22	Ground	Pin	Ground	Ground																																																																																																																																																						
C	23	Spare	Pin																																																																																																																																																								
C	24	Spare	Pin																																																																																																																																																								
EXCITE HD (12.x) and HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM HYPERTRONIC CONNECTOR PINOUTS (12.x – Forward Production & Upgrades) (14.x & 15.x – Upgrades) SHEET 9F OF 18																																																																																																																																																											

Signal HD, HDx & HDxt 3.0T (12.x, 14.x, 15.x & 16.x) Block Diagrams & Supplemental Schematics																																																											
REV 12																																																											
DIRECTION 5133320																																																											
HYPERTRONICS CONNECTOR (ROW C PINOUTS) DC, COIL ID, AND MULTICOIL TO SYSTEM INTERCONNECT																																																											
HYPERTRONICS CONNECTOR (ROW E PINOUTS) SYSTEM PREAMPLIFIER MULTICOIL TO SYSTEM INTERCONNECT																																																											
<table><tr><th colspan="2">Pin Out</th><th>Description</th><th>Style</th><th>PreAmp on Coil Connection Requirements</th><th>System PreAmp Connection Requirements</th></tr><tr><td>D</td><td>1</td><td>Coil Integrated Preamp RF1 In</td><td>Co–Axial</td><td>PreAmp on Coil RF in / No Bias</td><td></td></tr><tr><td>D</td><td>2</td><td>Coil Integrated Preamp RF2 In</td><td>Co–Axial</td><td>PreAmp on Coil RF in / No Bias</td><td></td></tr><tr><td>D</td><td>3</td><td>Coil Integrated Preamp RF3 In</td><td>Co–Axial</td><td>PreAmp on Coil RF in / No Bias</td><td></td></tr><tr><td>D</td><td>4</td><td>Coil Integrated Preamp RF4 In</td><td>Co–Axial</td><td>PreAmp on Coil RF in / No Bias</td><td></td></tr><tr><td>D</td><td>5</td><td>Coil Integrated Preamp RF5 In</td><td>Co–Axial</td><td>PreAmp on Coil RF in / No Bias</td><td></td></tr><tr><td>D</td><td>6</td><td>Coil Integrated Preamp RF6 In</td><td>Co–Axial</td><td>PreAmp on Coil RF in / No Bias</td><td></td></tr><tr><td>D</td><td>7</td><td>Coil Integrated Preamp RF7 In</td><td>Co–Axial</td><td>PreAmp on Coil RF in / No Bias</td><td></td></tr><tr><td>D</td><td>8</td><td>Coil Integrated Preamp RF8 In</td><td>Co–Axial</td><td>PreAmp on Coil RF in / No Bias</td><td></td></tr></table>						Pin Out		Description	Style	PreAmp on Coil Connection Requirements	System PreAmp Connection Requirements	D	1	Coil Integrated Preamp RF1 In	Co–Axial	PreAmp on Coil RF in / No Bias		D	2	Coil Integrated Preamp RF2 In	Co–Axial	PreAmp on Coil RF in / No Bias		D	3	Coil Integrated Preamp RF3 In	Co–Axial	PreAmp on Coil RF in / No Bias		D	4	Coil Integrated Preamp RF4 In	Co–Axial	PreAmp on Coil RF in / No Bias		D	5	Coil Integrated Preamp RF5 In	Co–Axial	PreAmp on Coil RF in / No Bias		D	6	Coil Integrated Preamp RF6 In	Co–Axial	PreAmp on Coil RF in / No Bias		D	7	Coil Integrated Preamp RF7 In	Co–Axial	PreAmp on Coil RF in / No Bias		D	8	Coil Integrated Preamp RF8 In	Co–Axial	PreAmp on Coil RF in / No Bias	
Pin Out		Description	Style	PreAmp on Coil Connection Requirements	System PreAmp Connection Requirements																																																						
D	1	Coil Integrated Preamp RF1 In	Co–Axial	PreAmp on Coil RF in / No Bias																																																							
D	2	Coil Integrated Preamp RF2 In	Co–Axial	PreAmp on Coil RF in / No Bias																																																							
D	3	Coil Integrated Preamp RF3 In	Co–Axial	PreAmp on Coil RF in / No Bias																																																							
D	4	Coil Integrated Preamp RF4 In	Co–Axial	PreAmp on Coil RF in / No Bias																																																							
D	5	Coil Integrated Preamp RF5 In	Co–Axial	PreAmp on Coil RF in / No Bias																																																							
D	6	Coil Integrated Preamp RF6 In	Co–Axial	PreAmp on Coil RF in / No Bias																																																							
D	7	Coil Integrated Preamp RF7 In	Co–Axial	PreAmp on Coil RF in / No Bias																																																							
D	8	Coil Integrated Preamp RF8 In	Co–Axial	PreAmp on Coil RF in / No Bias																																																							
HYPERTRONICS CONNECTOR (ROW E PINOUTS) SYSTEM PREAMPLIFIER MULTICOIL TO SYSTEM INTERCONNECT																																																											
<table><tr><th colspan="2">Pin Out</th><th>Description</th><th>Style</th><th>PreAmp on Coil Connection Requirements</th><th>System PreAmp Connection Requirements</th></tr><tr><td>E</td><td>1</td><td>RF1 In to System PreAmp</td><td>Co–Axial</td><td></td><td>Non–PreAmp on Coil RF in with / Bias</td></tr><tr><td>E</td><td>2</td><td>RF2 In to System PreAmp</td><td>Co–Axial</td><td></td><td>Non–PreAmp on Coil RF in with / Bias</td></tr><tr><td>E</td><td>3</td><td>RF3 In to System PreAmp</td><td>Co–Axial</td><td></td><td>Non–PreAmp on Coil RF in with / Bias</td></tr><tr><td>E</td><td>4</td><td>RF4 In to System PreAmp</td><td>Co–Axial</td><td></td><td>Non–PreAmp on Coil RF in with / Bias</td></tr><tr><td>E</td><td>5</td><td>RF5 In to System PreAmp</td><td>Co–Axial</td><td></td><td>Non–PreAmp on Coil RF in with / Bias</td></tr><tr><td>E</td><td>6</td><td>RF6 In to System PreAmp</td><td>Co–Axial</td><td></td><td>Non–PreAmp on Coil RF in with / Bias</td></tr><tr><td>E</td><td>7</td><td>RF7 In to System PreAmp</td><td>Co–Axial</td><td></td><td>Non–PreAmp on Coil RF in with / Bias</td></tr><tr><td>E</td><td>8</td><td>RF8 In to System PreAmp</td><td>Co–Axial</td><td></td><td>Non–PreAmp on Coil RF in with / Bias</td></tr></table>						Pin Out		Description	Style	PreAmp on Coil Connection Requirements	System PreAmp Connection Requirements	E	1	RF1 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias	E	2	RF2 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias	E	3	RF3 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias	E	4	RF4 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias	E	5	RF5 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias	E	6	RF6 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias	E	7	RF7 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias	E	8	RF8 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias
Pin Out		Description	Style	PreAmp on Coil Connection Requirements	System PreAmp Connection Requirements																																																						
E	1	RF1 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias																																																						
E	2	RF2 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias																																																						
E	3	RF3 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias																																																						
E	4	RF4 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias																																																						
E	5	RF5 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias																																																						
E	6	RF6 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias																																																						
E	7	RF7 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias																																																						
E	8	RF8 In to System PreAmp	Co–Axial		Non–PreAmp on Coil RF in with / Bias																																																						
MGD/RRF/RF																																																											
1–32																																																											

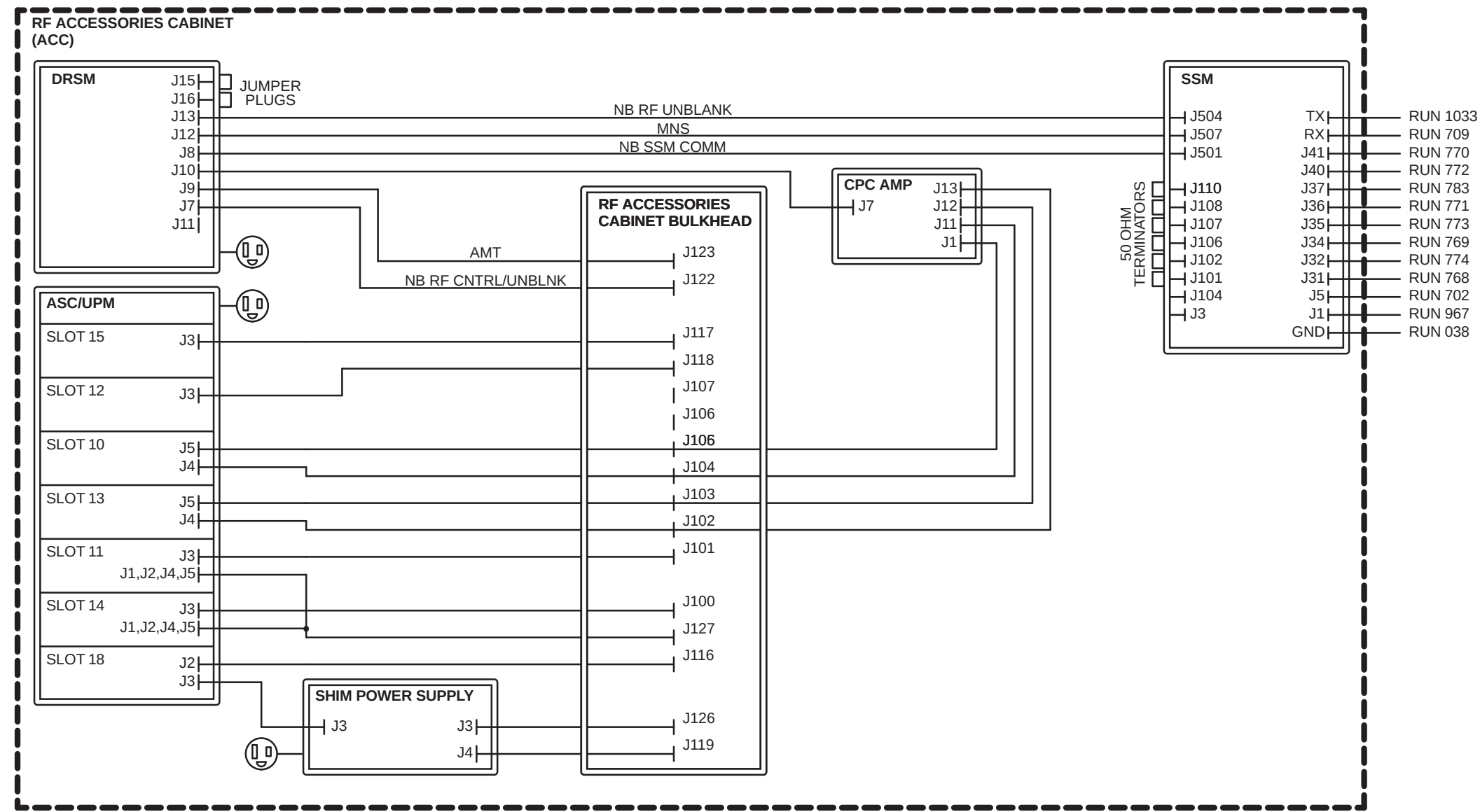

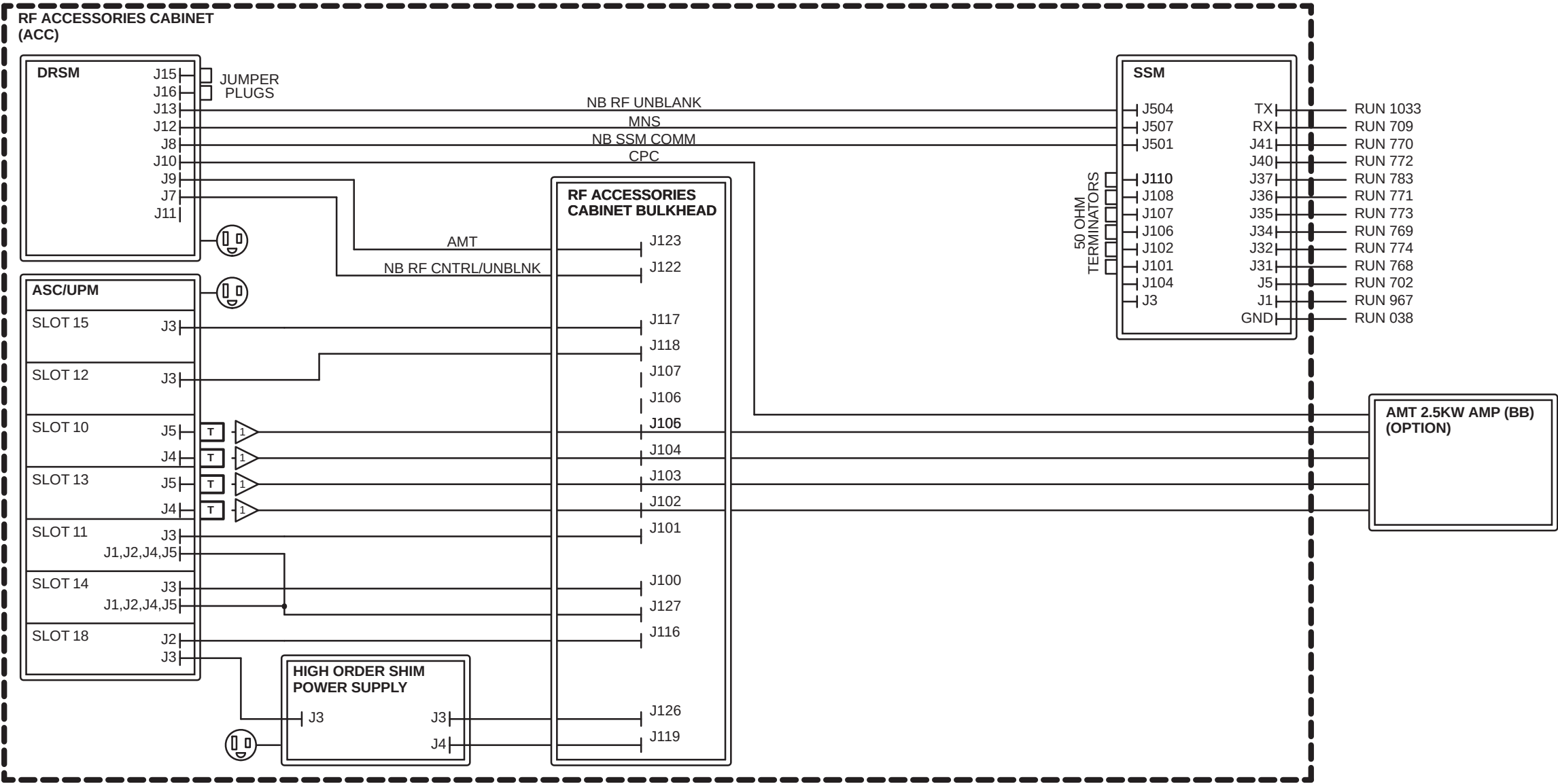

NOTES:
1 IF SITE HAS MNS, THESE CONNECTIONS WILL GO TO MNS AMPLIFIER CABINET. IF MNS NOT PART OF SYSTEM, THEN UPM BOARDS MUST HAVE 50 OHM TERMINATORS AT J4 AND J5 AT SLOTS 10 AND 13.

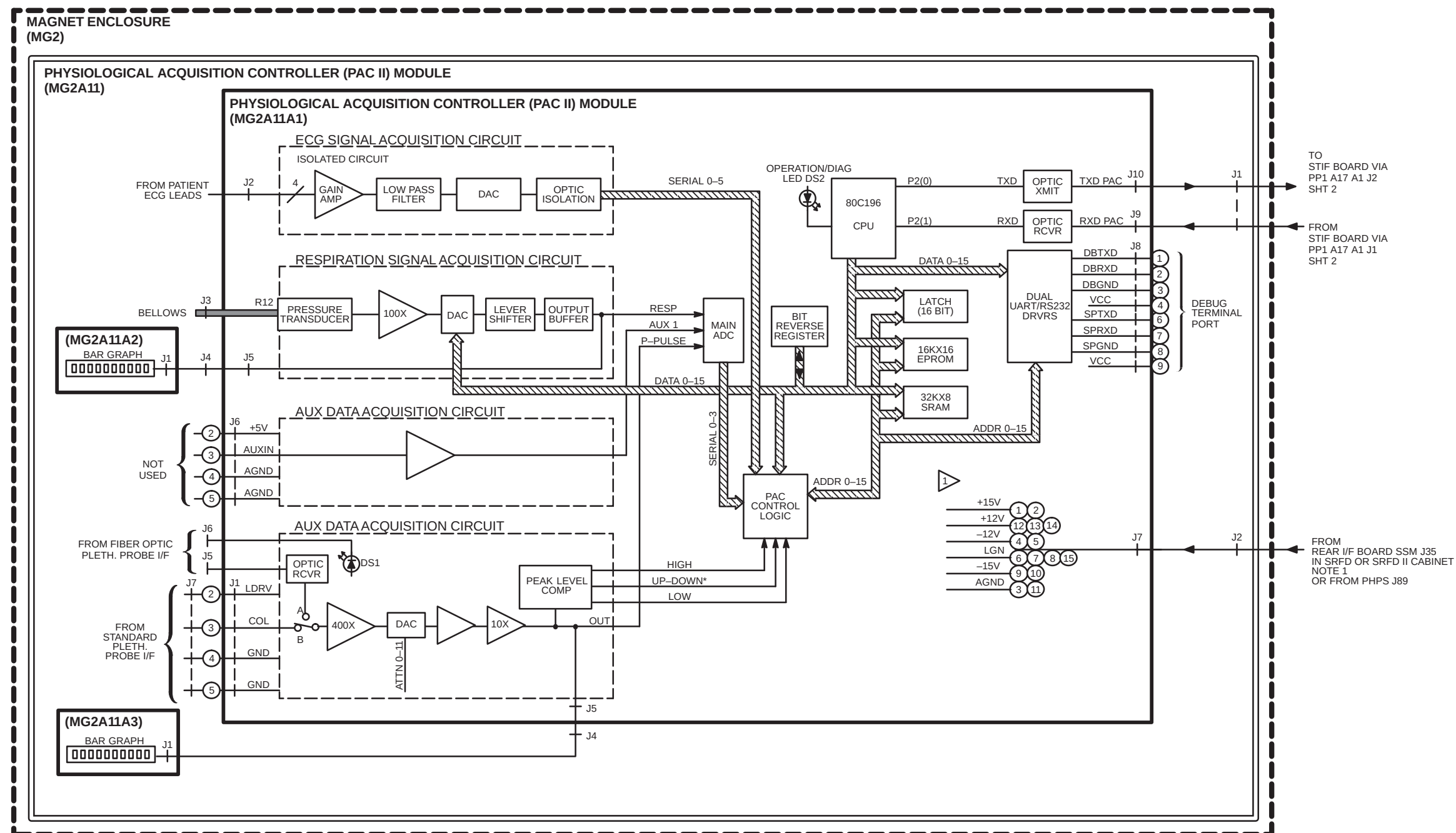

NOTES:



1 ▶ MULTIPLE WIRES/PINS WITH SAME SIGNAL ARE NOT SHOWN

REFER TO THE FOLLOWING FUNCTIONAL BLOCK DIAGRAM FOR ADDITIONAL INFORMATION:

NOTE:
1) PATIENT HANDLING

G3 & 3T94 EXCITE HD (12.x) and HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
PHYSIOLOGICAL ACQUISITION CONTROLLER II

SHEET 12 OF 18

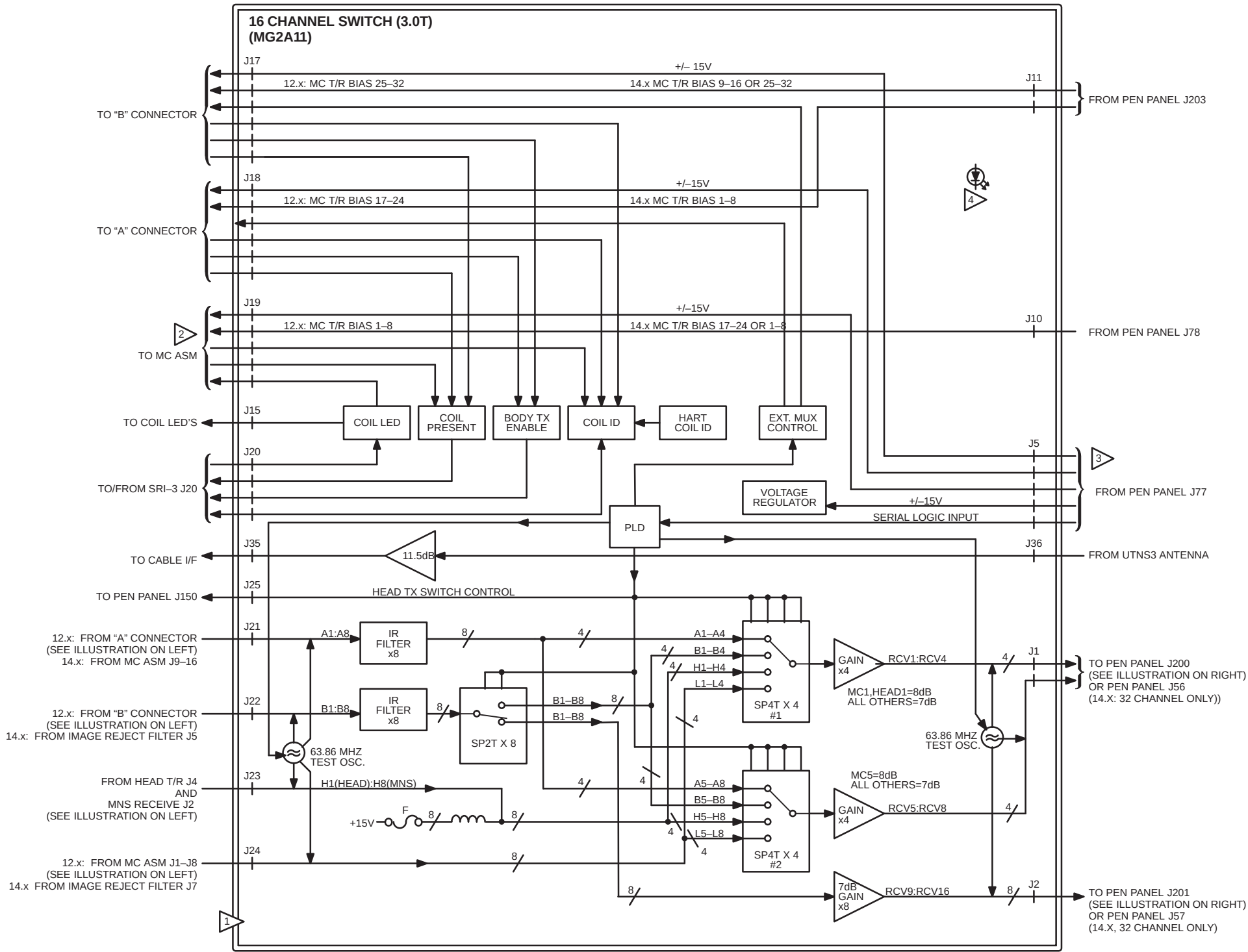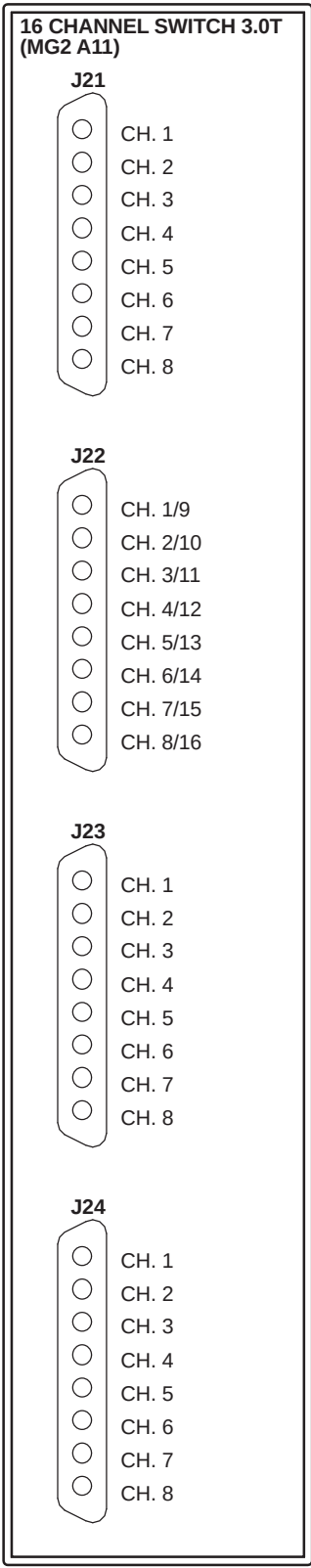

NOTES:

- Switch positions shown in 16 Channel Receive mode.
- More connector pinouts present on sheet 9A to 9C.

J5 PINS . SIGNALS

- 1 XPT_DATA_POS
- 2 XPT_STROBE_POS
- 3 XPT_OE_POS
- 4, 5 OPEN
- 6, 7, 8 ... +15V
- 9, 10 OPEN
- 11, 12 ... GND
- 13 OPEN
- 14 XPT_DATA_N
- 15 XPT_STROBE_N
- 16 XPT_OE_N_N
- 17, 18 ... OPEN
- 19, 20, 21 GND
- 22, 23 ... OPEN
- 24, 25 ... -15V

LED's

- COIL ID A GREEN
- COIL ID A RED
- COIL ID B GREEN
- COIL ID B RED
- COIL ID LEGACY GREEN
- COIL ID LEGACY RED
- COIL PRESENT A
- COIL PRESENT B
- COIL PRESENT 1
- COIL PRESENT 2
- TX ENABLE A
- TX ENABLE B
- DIAG OUT
- DIAG IN
- -15V POWER
- +15V POWER
- +24V POWER
- J21 ON CH. 1-4
- J22 ON CH. 1-4
- J23 ON CH. 1-4
- J24 ON CH. 1-4
- J21 ON CH. 5-8
- J22 ON CH. 5-8
- J23 ON CH. 5-8
- J24 ON CH. 5-8
- J22 ON CH. 9-12
- J22 ON CH. 13-16



16 CHANNEL SWITCH 1.5T (MG2 A11)

J1


- CH. 1
- CH. 2
- CH. 3
- CH. 4
- CH. 5
- CH. 6
- CH. 7
- CH. 8


J2

- CH. 9
- CH. 10
- CH. 11
- CH. 12
- CH. 13
- CH. 14
- CH. 15
- CH. 16

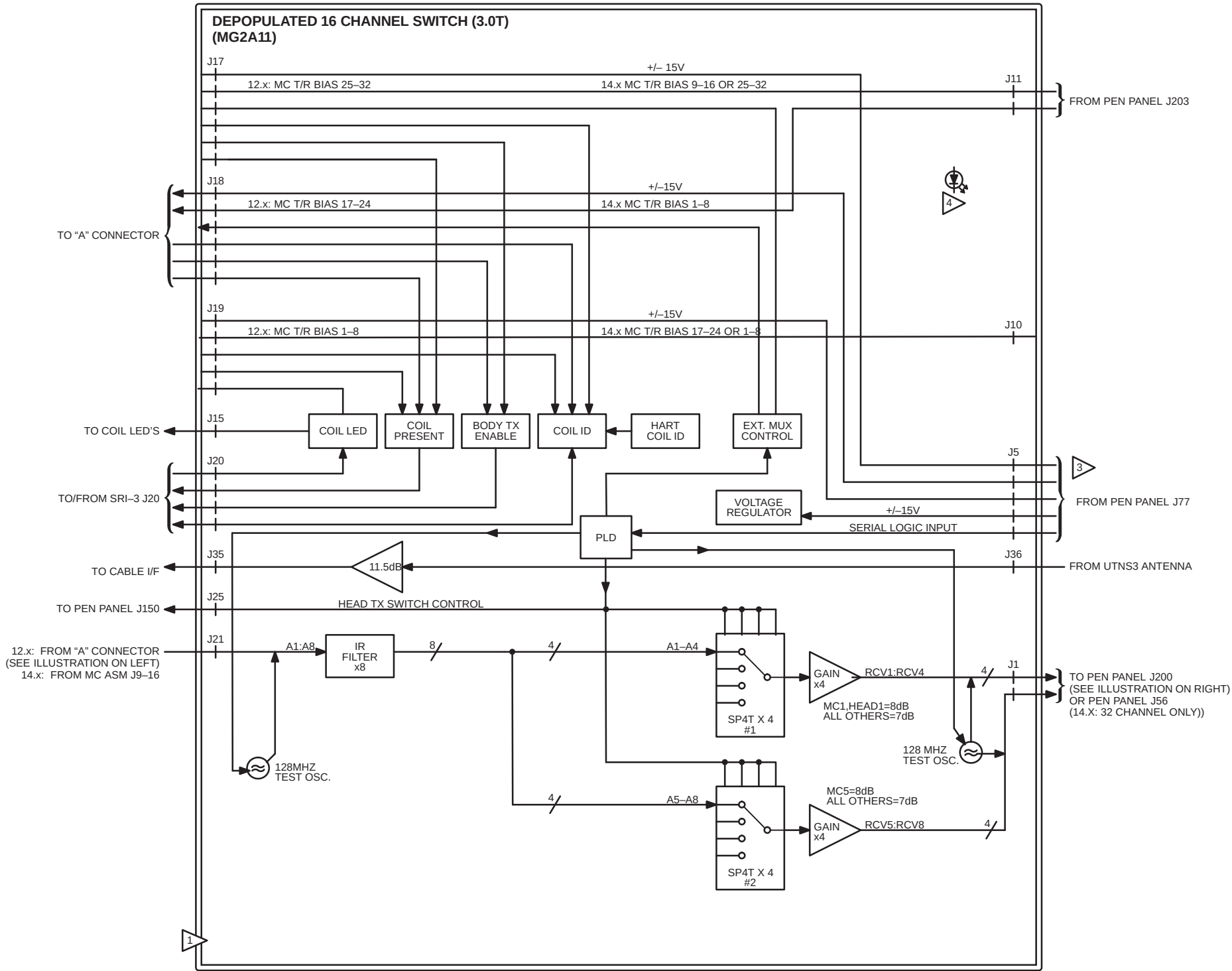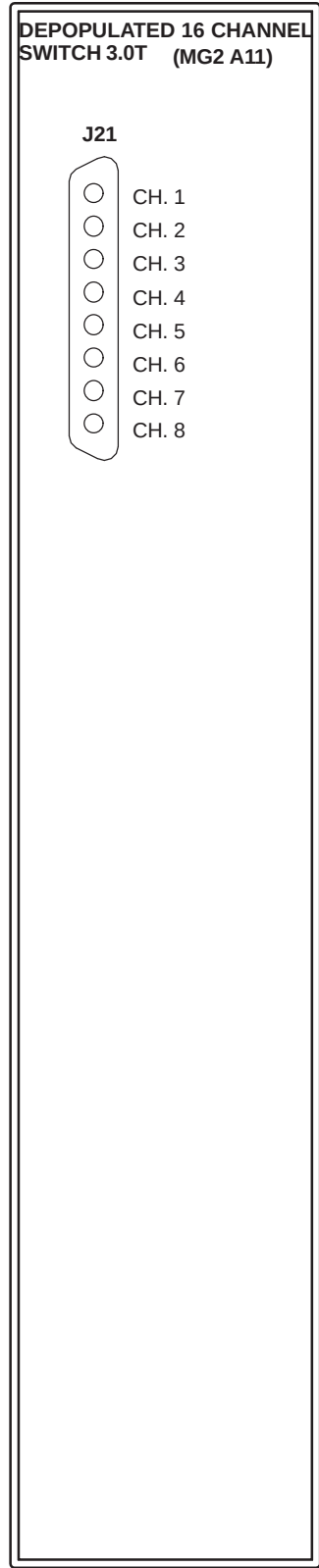


- NOTES:
- 1 Switch positions shown in 8 Channel Receive mode.
- 2 More connector pinouts present on sheet 9A.
- 3 J5 PINS . SIGNALS
- | | | |
|------------|-------|----------------|
| 1 | | XPT_DATA_POS |
| 2 | | XPT_STROBE_POS |
| 3 | | XPT_OE_POS |
| 4, 5 | | OPEN |
| 6, 7, 8 | ... | +15V |
| 9, 10 | ... | OPEN |
| 11, 12 | ... | GND |
| 13 | | OPEN |
| 14 | | XPT_DATA_N |
| 15 | | XPT_STROBE_N |
| 16 | | XPT_OE_N |
| 17, 18 | ... | OPEN |
| 19, 20, 21 | GND | |
| 22, 23 | ... | OPEN |
| 24, 25 | ... | -15V |
- 4 LED'S
- | | |
|-------------------|------------------|
| • COIL ID A GREEN | • -15V POWER |
| • COIL ID A RED | • +15V POWER |
| • COIL PRESENT A | • +24V POWER |
| • TX ENABLE A | • J21 ON CH. 1-4 |

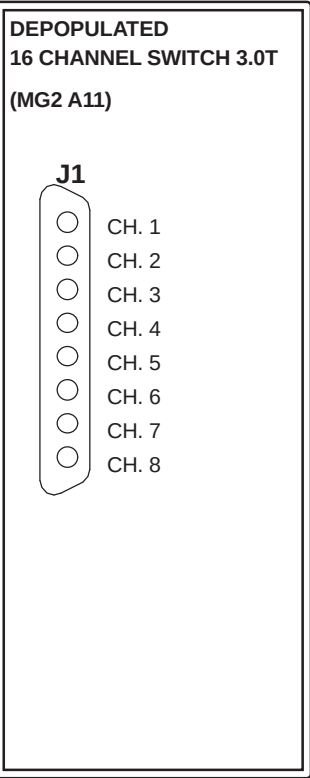

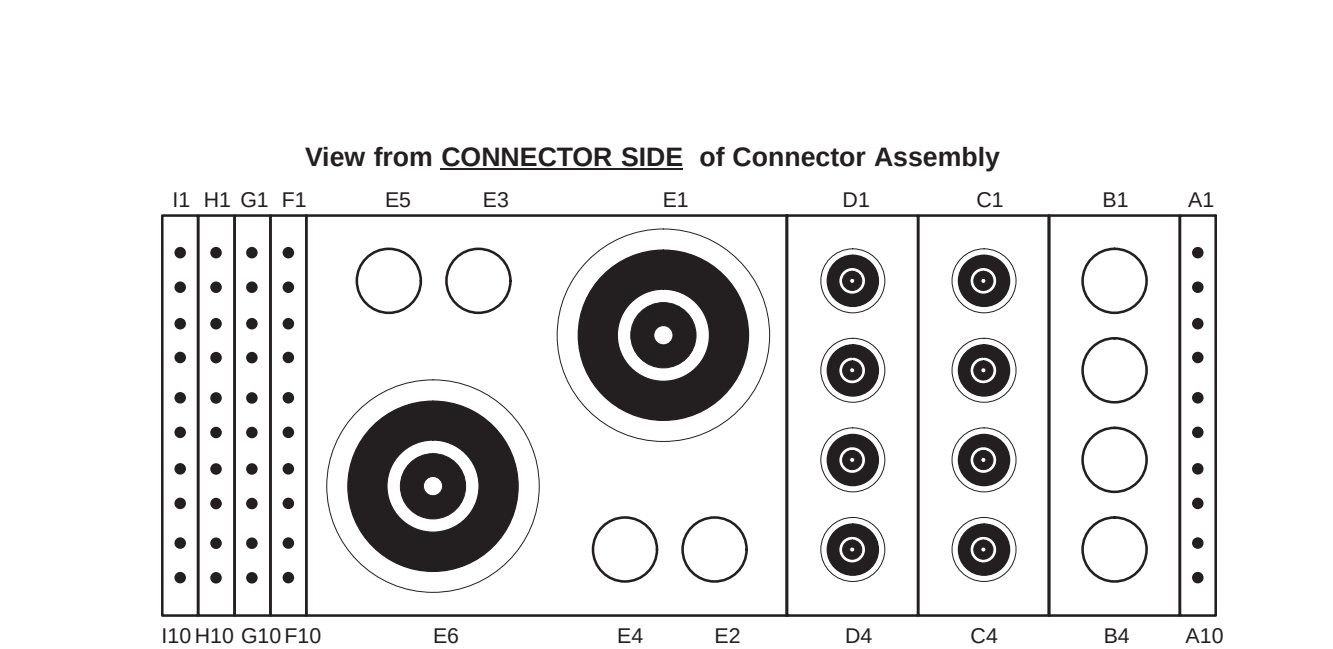


Contact Function Allocation			
Pin	Twisted with	Description	Comments
A1	A2	coil present	
A2	A1	coil present return	
A3	A4	multiplexerSelect0	RS-422+
A4	A3	multiplexerSelect0 return	RS-422-
A5	A6	multiplexerSelect1	RS-422+
A6	A5	multiplexerSelect1 return	RS-422-
A7	A8	multiplexerSelect2	RS-422+
A8	A7	multiplexerSelect2 return	RS-422-
A9	A10	multiplexerSelect3	RS-422+
A10	A9	multiplexerSelect3 return	RS-422-
B1	n/a	Key blocked in receptacle E	
B2	n/a	Key blocked in receptacle F	
B3	n/a	Key blocked in receptacle G	
B4	n/a	Key blocked in receptacle H	
C1	n/a	MC-1A,B or MC-9B	Preamp out 1, no bias, DC blocked
C2	n/a	MC-2A,B or MC-10B	Preamp out 2, no bias, DC blocked
C3	n/a	MC-3A,B or MC-11B	Preamp out 3, no bias, DC blocked
C4	n/a	MC-4A,B or MC-12B	Preamp out 4, no bias, DC blocked
D1	n/a	MC-5A,B or MC-13B	Preamp out 5, no bias, DC blocked
D2	n/a	MC-6A,B or MC-14B	Preamp out 6, no bias, DC blocked
D3	n/a	MC-7A,B or MC-15B	Preamp out 7, no bias, DC blocked
D4	n/a	MC-8A,B or MC-16B	Preamp out 8, no bias, DC blocked

Contact Function Allocation			
Pin	Twisted with	Description	Comments
E1	n/a	Transmit signal up to 2kw (1.5T) or 4kw (3.0T) PEP.	“A” Connector only. Not present in “B” Connector
E2	n/a	Key blocked in receptacle C	
E3	n/a	Key blocked in receptacle A	Blocked for “A”
E4	n/a	Key blocked in receptacle B	Blocked for “B”
E5	n/a	Key blocked in receptacle D	
E6	n/a	Reflected power to 50 ohm load 2.5kw PEP, 50w RMS.	“A” Connector only. Not present in “B” Connector
F1	G1	MC bias 1	MC T/R Bias Line 17 for “A”, 25 for “B”
F2	G2	MC bias 2	MC T/R Bias Line 18 for “A”, 26 for “B”
F3	G3	MC bias 3	MC T/R Bias Line 19 for “A”, 27 for “B”
F4	G4	MC bias 4	MC T/R Bias Line 20 for “A”, 28 for “B”
F5	G5	MC bias 5	MC T/R Bias Line 21 for “A”, 29 for “B”
F6	G6	MC bias 6	MC T/R Bias Line 22 for “A”, 30 for “B”
F7	G7	MC bias 7	MC T/R Bias Line 23 for “A”, 31 for “B”
F8	G8	MC bias 8	MC T/R Bias Line 24 for “A”, 32 for “B”
F9	F10	body transmit enable	
F10	F9	body transmit enable return	
G1	F1	MC bias return 1	MC T/R Bias Line 17 for “A”, 25 for “B”
G2	F2	MC bias return 2	MC T/R Bias Line 18 for “A”, 26 for “B”
G3	F3	MC bias return 3	MC T/R Bias Line 19 for “A”, 27 for “B”
G4	F4	MC bias return 4	MC T/R Bias Line 20 for “A”, 28 for “B”
G5	F5	MC bias return 5	MC T/R Bias Line 21 for “A”, 29 for “B”
G6	F6	MC bias return 6	MC T/R Bias Line 22 for “A”, 30 for “B”
G7	F7	MC bias return 7	MC T/R Bias Line 23 for “A”, 31 for “B”
G8	F8	MC bias return 8	MC T/R Bias Line 24 for “A”, 32 for “B”
G9	n/a	Reserved	
G10	n/a	Reserved	
H1	n/a	Reserved	
H2	n/a	Reserved	
H3	n/a	Reserved	
H4	n/a	Reserved	
H5	n/a	Reserved	
H6	n/a	Reserved	
H7	n/a	Reserved	
H8	n/a	Reserved	
H9	n/a	Reserved	
H10	n/a	Reserved	
I1	n/a	Reserved	
I2	I5	+15Vdc (auxiliary & preamp power)	
I3	I6	+10Vdc (auxiliary & preamp power)	
I4	n/a	-15Vdc (auxiliary power)	
I5	I2	Return for aux power	
I6	I3	Return for aux power	
I7	I8	coil ID (Dallas OneWire protocol)	All coils MUST have coil ID.
I8	I7	coil ID return (Dallas OneWire protocol)	All coils MUST have coil ID.
I9	I10	coil present'	
I10	I9	coil present return'	

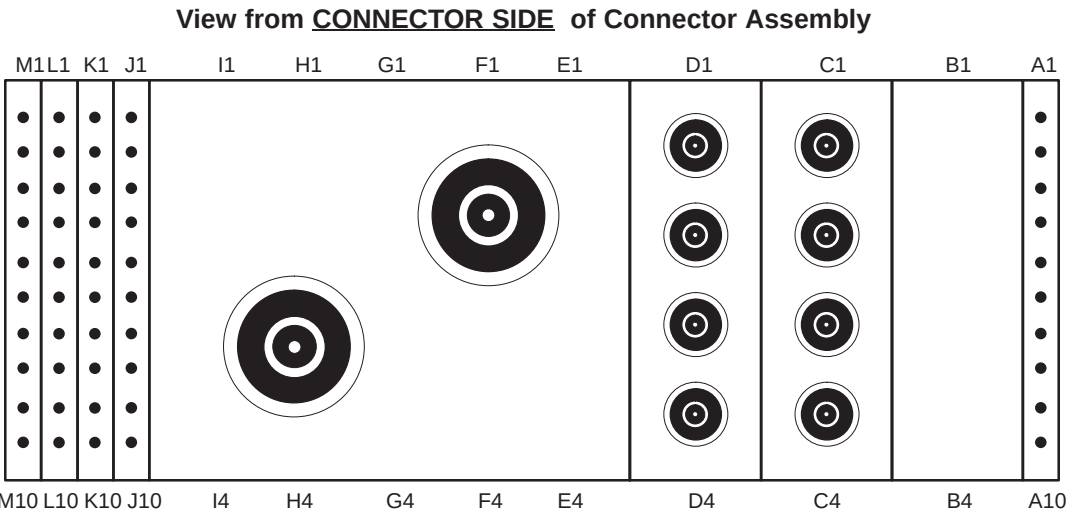


Contact Function Allocation			
Pin	Twisted with	Description	Comments
A1	A2	coil present	
A2	A1	coil present return	
A3	A4	multiplexerSelect0	RS-422+
A4	A3	multiplexerSelect0 return	RS-422-
A5	A6	multiplexerSelect1	RS-422+
A6	A5	multiplexerSelect1 return	RS-422-
A7	A8	multiplexerSelect2	RS-422+
A8	A7	multiplexerSelect2 return	RS-422-
A9	A10	multiplexerSelect3	RS-422+
A10	A9	multiplexerSelect3 return	RS-422-
B1	n/a	Key blocked in receptacle E	
B2	n/a	Key blocked in receptacle F	
B3	n/a	Key blocked in receptacle G	
B4	n/a	Key blocked in receptacle H	
C1	n/a	MC-A1	Preamp out 1, no bias, DC blocked
C2	n/a	MC-A2	Preamp out 2, no bias, DC blocked
C3	n/a	MC-A3	Preamp out 3, no bias, DC blocked
C4	n/a	MC-A4	Preamp out 4, no bias, DC blocked
D1	n/a	MC-A5	Preamp out 5, no bias, DC blocked
D2	n/a	MC-A6	Preamp out 6, no bias, DC blocked
D3	n/a	MC-A7	Preamp out 7, no bias, DC blocked
D4	n/a	MC-A8	Preamp out 8, no bias, DC blocked

Contact Function Allocation			
Pin	Twisted with	Description	Comments
F2	n/a	Transmit Signal up to 8kw PEP	
H3	n/a	Reflected power to 50 ohm load 2.5kw PEP, 50w RMS	"A" connector only
J1	K1	MC bias 1	
J2	K2	MC bias 2	
J3	K3	MC bias 3	
J4	K4	MC bias 4	
J5	K5	MC bias 5	
J6	K6	MC bias 6	
J7	K7	MC bias 7	
J8	K8	MC bias 8	
J9	J10	body transmit enable	
J10	J9	body transmit enable return	
K1	J1	MC bias return 1	
K2	J2	MC bias return 2	
K3	J3	MC bias return 3	
K4	J4	MC bias return 4	
K5	J5	MC bias return 5	
K6	J6	MC bias return 6	
K7	J7	MC bias return 7	
K8	J8	MC bias return 8	
K9	K10	System preamp enable return	
K10	K9	System preamp enable	
L1	n/a	Reserved	
L2	n/a	Reserved	
L3	n/a	Reserved	
L4	n/a	Reserved	
L5	n/a	Reserved	
L6	n/a	Reserved	
L7	n/a	Reserved	
L8	n/a	Reserved	
L9	n/a	Reserved	
L10	n/a	Reserved	
M1	n/a	Reserved	
M2	M5	+15Vdc (auxiliary & preamp power)	
M3	M6	+10Vdc (auxiliary & preamp power)	
M4	n/a	-15Vdc (auxiliary power)	
M5	M2	Return for aux power	
M6	M3	Return for aux power	
M7	M8	coil ID (Dallas OneWire protocol)	All coils MUST have coil ID.
M8	M7	coil ID return (Dallas OneWire protocol)	All coils MUST have coil ID.
M9	M10	coil present'	
M10	M9	coil present return'	

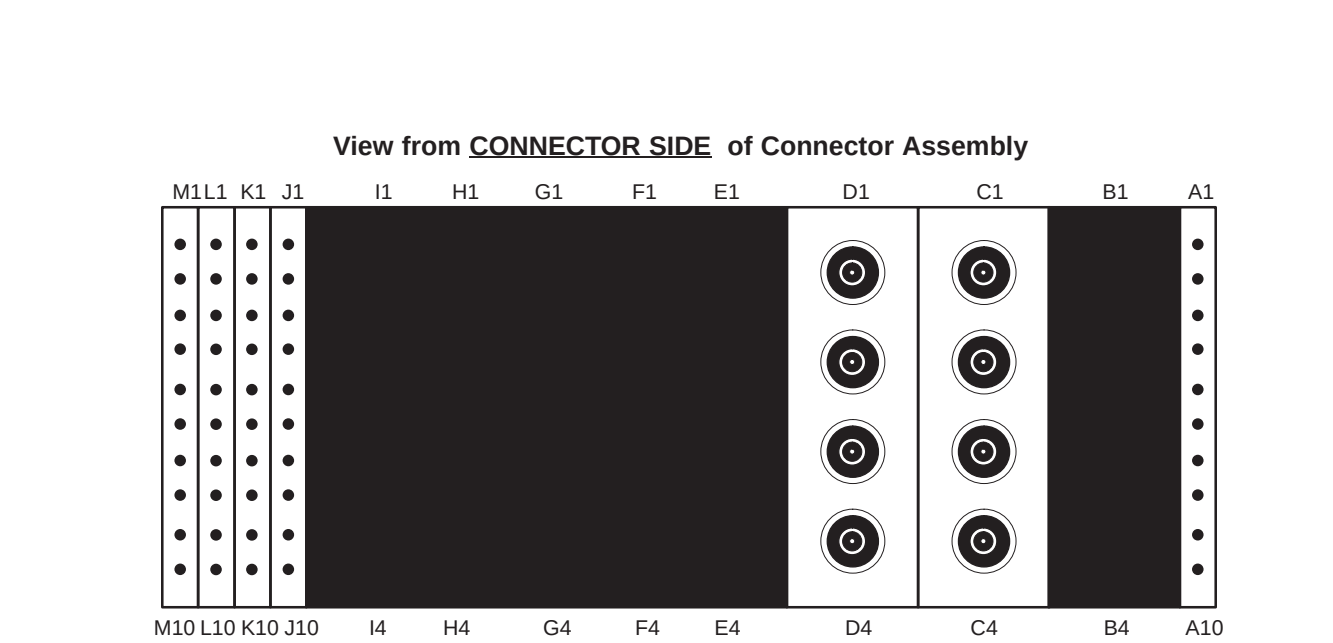


Contact Function Allocation			
Pin	Twisted with	Description	Comments
A1	A2	coil present	
A2	A1	coil present return	
A3	A4	multiplexerSelect0	RS-422+
A4	A3	multiplexerSelect0 return	RS-422-
A5	A6	multiplexerSelect1	RS-422+
A6	A5	multiplexerSelect1 return	RS-422-
A7	A8	multiplexerSelect2	RS-422+
A8	A7	multiplexerSelect2 return	RS-422-
A9	A10	multiplexerSelect3	RS-422+
A10	A9	multiplexerSelect3 return	RS-422-
B1	n/a		
B2	n/a		
B3	n/a		
B4	n/a		
C1	n/a	MC-B1	Preamp out 1, MC bias 1 (9)
C2	n/a	MC-B2	Preamp out 2, MC bias 2 (10)
C3	n/a	MC-B3	Preamp out 3, MC bias 3 (11)
C4	n/a	MC-B4	Preamp out 4, MC bias 4 (12)
D1	n/a	MC-B5	Preamp out 5, MC bias 5 (13)
D2	n/a	MC-B6	Preamp out 6, MC bias 6 (14)
D3	n/a	MC-B7	Preamp out 7, MC bias 7 (15)
D4	n/a	MC-B8	Preamp out 8, MC bias 8 (16)

Contact Function Allocation			
Pin	Twisted with	Description	Comments
F2	n/a		
J1	K1	MC bias 1 (9)	
J2	K2	MC bias 2 (10)	
J3	K3	MC bias 3 (11)	
J4	K4	MC bias 4 (12)	
J5	K5	MC bias 5 (13)	
J6	K6	MC bias 6 (14)	
J7	K7	MC bias 7 (15)	
J8	K8	MC bias 8 (16)	
J9	J10	body transmit enable	
J10	J9	body transmit enable return	
K1	J1	MC bias return 1 (9)	
K2	J2	MC bias return 2 (10)	
K3	J3	MC bias return 3 (11)	
K4	J4	MC bias return 4 (12)	
K5	J5	MC bias return 5 (13)	
K6	J6	MC bias return 6 (14)	
K7	J7	MC bias return 7 (15)	
K8	J8	MC bias return 8 (16)	
K9	K10	System preamp enable return	
K10	K9	System preamp enable	
L1	n/a	Reserved	
L2	n/a	Reserved	
L3	n/a	Reserved	
L4	n/a	Reserved	
L5	n/a	Reserved	
L6	n/a	Reserved	
L7	n/a	Reserved	
L8	n/a	Reserved	
L9	n/a	Reserved	
L10	n/a	Reserved	
M1	n/a	Reserved	
M2	M5	+15Vdc (auxiliary & preamp power)	
M3	M6	+10Vdc (auxiliary & preamp power)	
M4	n/a	-15Vdc (auxiliary power)	
M5	M2	Return for aux power	
M6	M3	Return for aux power	
M7	M8	coil ID (Dallas OneWire protocol)	All coils MUST have coil ID.
M8	M7	coil ID return (Dallas OneWire protocol)	All coils MUST have coil ID.
M9	M10	coil present'	
M10	M9	coil present return'	

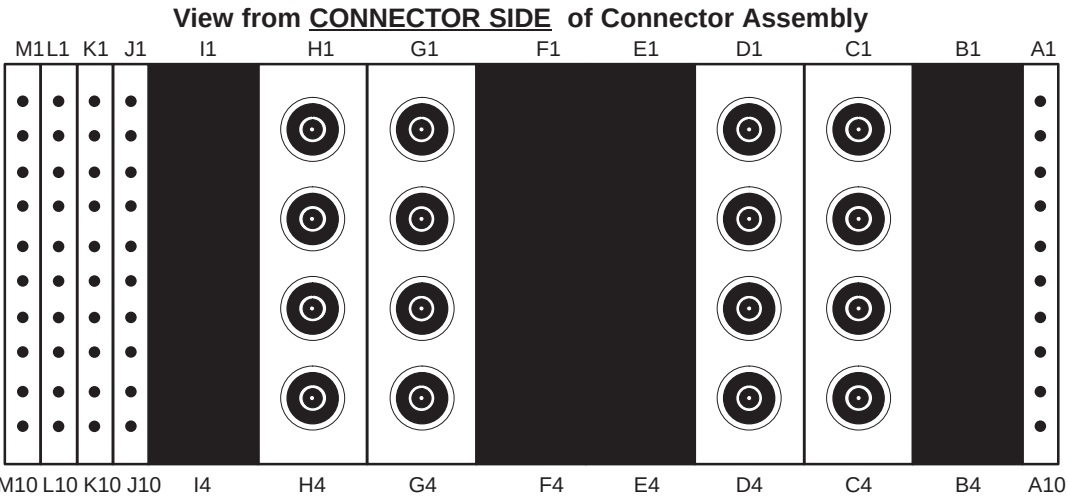



Contact Function Allocation			
Pin	Twisted with	Description	Comments
A1	A2	coil present	
A2	A1	coil present return	
A3	A4	multiplexerSelect0	RS-422+
A4	A3	multiplexerSelect0 return	RS-422-
A5	A6	multiplexerSelect1	RS-422+
A6	A5	multiplexerSelect1 return	RS-422-
A7	A8	multiplexerSelect2	RS-422+
A8	A7	multiplexerSelect2 return	RS-422-
A9	A10	multiplexerSelect3	RS-422+
A10	A9	multiplexerSelect3 return	RS-422-
C1	n/a	MC-B1 or MC-C1, MC bias 1B (9), 1C (25)	Preamp out 1, MC bias 1B (9), 1C (25)
C2	n/a	MC-B2 or MC-C2, MC bias 2B (10), 2C (26)	Preamp out 2, MC bias 2B (10), 2C (26)
C3	n/a	MC-B3 or MC-C3, MC bias 3B (11), 3C (27)	Preamp out 3, MC bias 3B (11), 3C (27)
C4	n/a	MC-B4 or MC-C4, MC bias 4B (12), 4C (28)	Preamp out 4, MC bias 4B (12), 4C (28)
D1	n/a	MC-B5 or MC-C5, MC bias 5B (13), 5C (29)	Preamp out 5, MC bias 5B (13), 5C (29)
D2	n/a	MC-B6 or MC-C6, MC bias 6B (14), 6C (30)	Preamp out 6, MC bias 6B (14), 6C (30)
D3	n/a	MC-B7 or MC-C7, MC bias 7B (15), 7C (31)	Preamp out 7, MC bias 7B (15), 7C (31)
D4	n/a	MC-B8 or MC-C8, MC bias 8B (16), 8C (32)	Preamp out 8, MC bias 8B (16), 8C (32)
G1	n/a	MC-B9 or MC-C9, MC bias 9B (25), 9C (1)	Preamp out 1, MC bias 9B (25), 9C (1)
G2	n/a	MC-B10 or MC-102, MC bias 10B (26), 10C (2)	Preamp out 2, MC bias 10B (26), 10C (2)
G3	n/a	MC-B11 or MC-C11, MC bias 11B (27), 11C (3)	Preamp out 3, MC bias 11B (27), 11C (3)
G4	n/a	MC-B12 or MC-C12, MC bias 12B (28), 12C (4)	Preamp out 4, MC bias 12B (28), 12C (4)
H1	n/a	MC-B13 or MC-C13, MC bias 13B (29), 13C (5)	Preamp out 5, MC bias 13B (29), 13C (5)
H2	n/a	MC-B14 or MC-C14, MC bias 14B (30), 14C (6)	Preamp out 6, MC bias 14B (30), 14C (6)
H3	n/a	MC-B15 or MC-C15, MC bias 15B (31), 15C (7)	Preamp out 7, MC bias 15B (31), 15C (7)
H4	n/a	MC-B16 or MC-C16, MC bias 16B (32), 16C (8)	Preamp out 8, MC bias 16B (32), 16C (8)

Contact Function Allocation			
Pin	Twisted with	Description	Comments
J1	K1	MC bias 1 (9)-B / MC bias 1 (25)-C	
J2	K2	MC bias 2 (10)-B / MC bias 2 (26)-C	
J3	K3	MC bias 3 (11)-B / MC bias 3 (27)-C	
J4	K4	MC bias 4 (12)-B / MC bias 4 (28)-C	
J5	K5	MC bias 5 (13)-B / MC bias 5 (29)-C	
J6	K6	MC bias 6 (14)-B / MC bias 6 (30)-C	
J7	K7	MC bias 7 (15)-B / MC bias 7 (31)-C	
J8	K8	MC bias 8 (16)-B / MC bias 8 (32)-C	
J9	J10	body transmit enable	
J10	J9	body transmit enable return	
K1	J1	MC bias return 1 (25)-B / MC bias return 1 (1)-C	
K2	J2	MC bias return 2 (26)-B / MC bias return 2 (2)-C	
K3	J3	MC bias return 3 (27)-B / MC bias return 3 (3)-C	
K4	J4	MC bias return 4 (28)-B / MC bias return 4 (4)-C	
K5	J5	MC bias return 5 (29)-B / MC bias return 5 (5)-C	
K6	J6	MC bias return 6 (30)-B / MC bias return 6 (6)-C	
K7	J7	MC bias return 7 (31)-B / MC bias return 7 (7)-C	
K8	J8	MC bias return 8 (32)-B / MC bias return 8 (8)-C	
K9	K10	System preamp enable return	
K10	K9	System preamp enable	
L1	n/a	Reserved	
L2	n/a	Reserved	
L3	n/a	Reserved	
L4	n/a	Reserved	
L5	n/a	Reserved	
L6	n/a	Reserved	
L7	n/a	Reserved	
L8	n/a	Reserved	
L9	n/a	Reserved	
L10	n/a	Reserved	
M1	n/a	Reserved	
M2	M5	+15Vdc (auxiliary & preamp power)	
M3	M6	+10Vdc (auxiliary & preamp power)	
M4	n/a	-15Vdc (auxiliary power)	
M5	M2	Return for aux power	
M6	M3	Return for aux power	
M7	M8	coil ID (Dallas OneWire protocol)	All coils MUST have coil ID.
M8	M7	coil ID return (Dallas OneWire protocol)	All coils MUST have coil ID.
M9	M10	coil present'	
M10	M9	coil present return'	

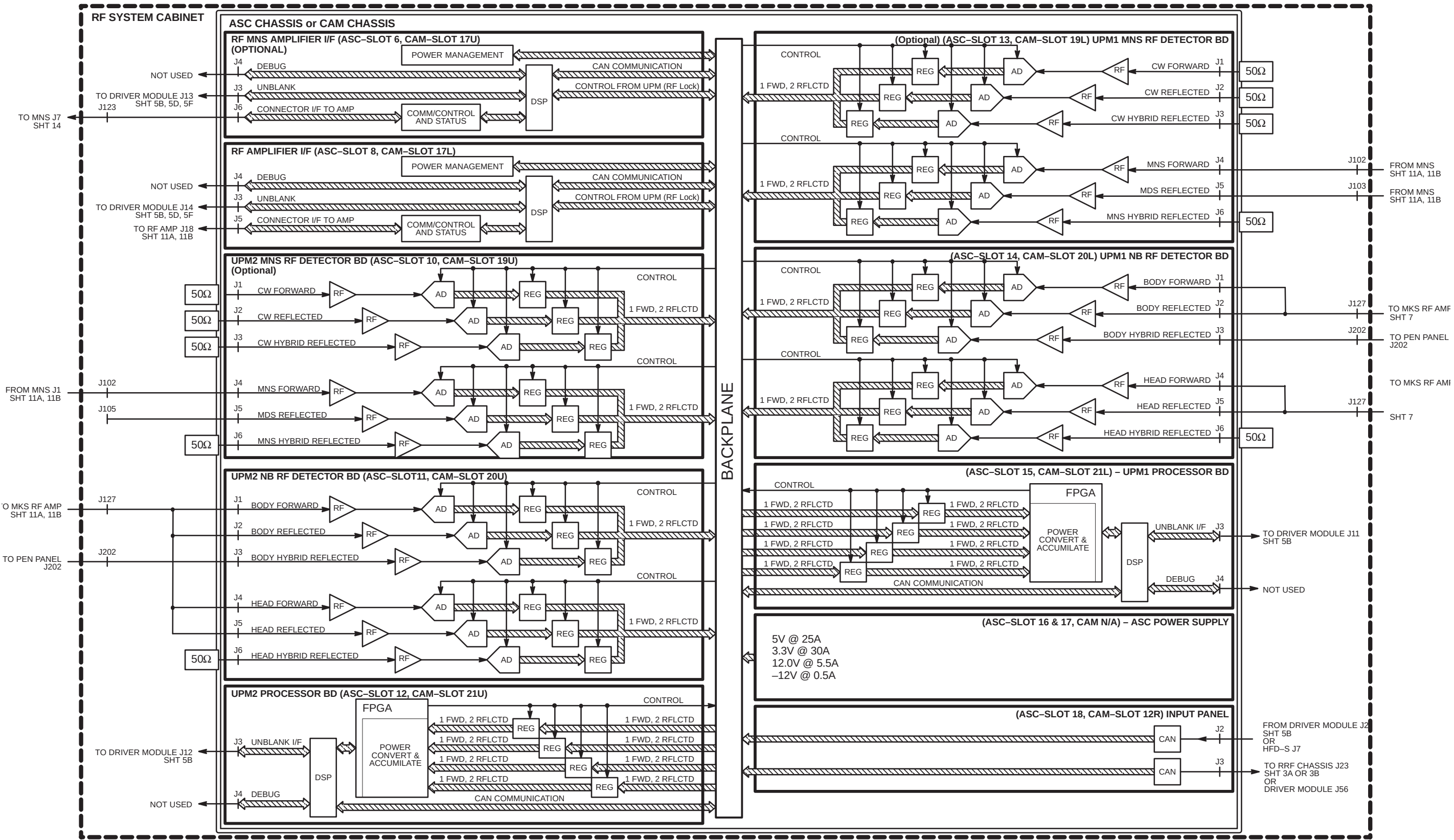


G3 & 3T94 EXCITE (12.x) and HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
ASC CHASSIS
SHEET 15 OF 18

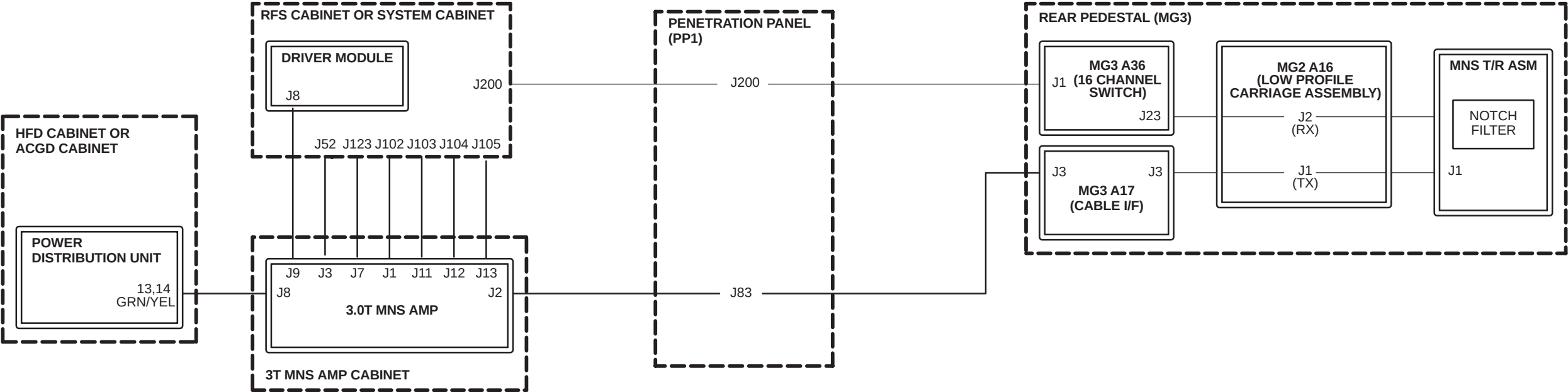

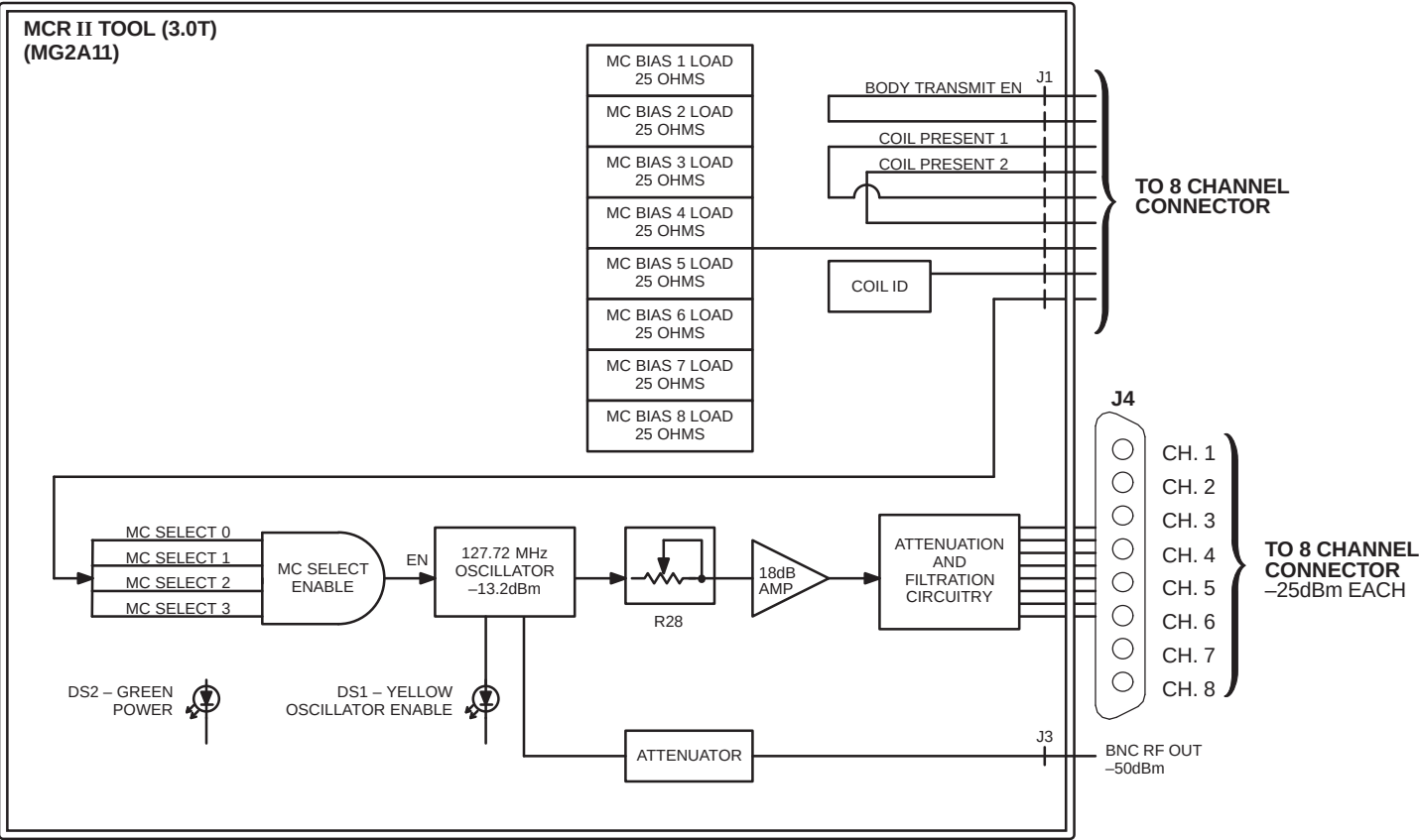


J1	
PINS	SIGNAL
1	MC Bias 1
2	MC Bias 2
3	MC Bias 3
4	MC Bias 4
5	MC Bias 5
6	MC Bias 6
7	MC Bias 7
8	MC Bias 8
9	MC Select 0
10	MC Select 1
11	MC Select 2
12	MC Select 3
13	Body Transmit Enable
14	Coil ID
15	Coil Present 15
16	Coil Present 2
17	+15VDC (Auxiliary Power)
18	+10VDC (Auxiliary Power)
19	-15VDC (Auxiliary Power)
20	MC Bias Return 1
21	MC Bias Return 2
22	MC Bias Return 3
23	MC Bias Return 4
24	MC Bias Return 5
25	MC Bias Return 6
26	MC Bias Return 7
27	MC Bias Return 8
28	MC Select 0 Return
29	MC Select 1 Return
30	MC Select 2 Return
31	MC Select 3 Return
32	Body Transmit Enable Return
33	Coil ID Return
34	Coil Present Return 1
35	Coil Present Return 2
36	Auxiliary Power Return
37	Auxiliary Power Return

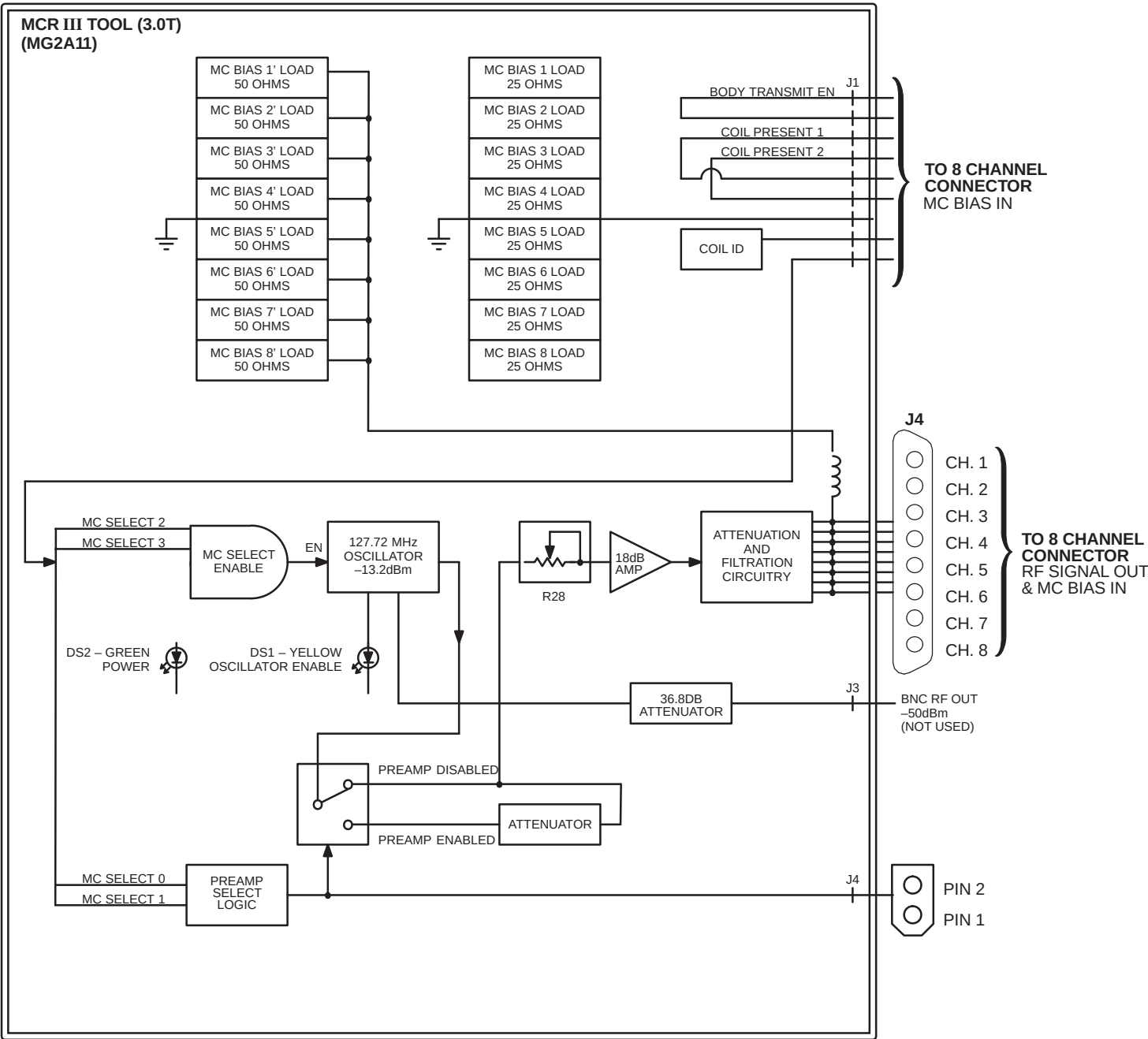

J1	
PINS	SIGNAL
1	MC Bias 1
2	MC Bias 2
3	MC Bias 3
4	MC Bias 4
5	MC Bias 5
6	MC Bias 6
7	MC Bias 7
8	MC Bias 8
9	MC Select 0
10	MC Select 1
11	MC Select 2
12	MC Select 3
13	Body Transmit Enable
14	Coil ID
15	Coil Present 15
16	Coil Present 2
17	+15VDC (Auxiliary Power)
18	+10VDC (Auxiliary Power)
19	-15VDC (Auxiliary Power)
20	MC Bias Return 1
21	MC Bias Return 2
22	MC Bias Return 3
23	MC Bias Return 4
24	MC Bias Return 5
25	MC Bias Return 6
26	MC Bias Return 7
27	MC Bias Return 8
28	MC Select 0 Return
29	MC Select 1 Return
30	MC Select 2 Return
31	MC Select 3 Return
32	Body Transmit Enable Return
33	Coil ID Return
34	Coil Present Return 1
35	Coil Present Return 2
36	Auxiliary Power Return
37	Auxiliary Power Return

J4	
PINS	SIGNAL
1	Preamp Enable
2	Preamp Enable RTN

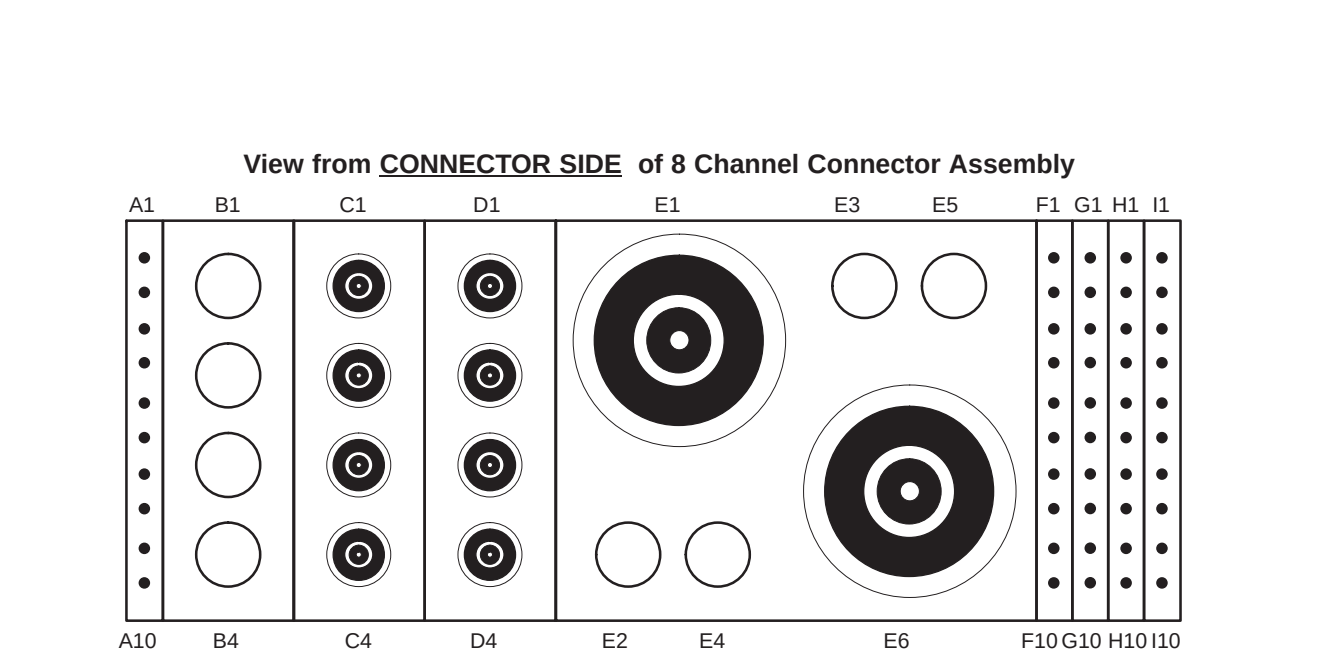


Contact Function Allocation			
Pin	Twisted with	Description	Comments
A1	A2	coil present	
A2	A1	coil present return	
A3	A4	multiplexerSelect0	RS-422+
A4	A3	multiplexerSelect0 return	RS-422-
A5	A6	multiplexerSelect1	RS-422+
A6	A5	multiplexerSelect1 return	RS-422-
A7	A8	multiplexerSelect2	RS-422+
A8	A7	multiplexerSelect2 return	RS-422-
A9	A10	multiplexerSelect3	RS-422+
A10	A9	multiplexerSelect3 return	RS-422-
B1	n/a		Not used
B2	n/a		Not used
B3	n/a		Not used
B4	n/a		Not used
C1	n/a	MC-1A,B or MC-9B	Preamp out 1, no bias, DC blocked
C2	n/a	MC-2A,B or MC-10B	Preamp out 2, no bias, DC blocked
C3	n/a	MC-3A,B or MC-11B	Preamp out 3, no bias, DC blocked
C4	n/a	MC-4A,B or MC-12B	Preamp out 4, no bias, DC blocked
D1	n/a	MC-5A,B or MC-13B	Preamp out 5, no bias, DC blocked
D2	n/a	MC-6A,B or MC-14B	Preamp out 6, no bias, DC blocked
D3	n/a	MC-7A,B or MC-15B	Preamp out 7, no bias, DC blocked
D4	n/a	MC-8A,B or MC-16B	Preamp out 8, no bias, DC blocked

G3 & 3T94 EXCITE HD (12.x) and HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
MCR II Tool – 8 Channel Connector
(12.x – Forward Production & Upgrades)
(14.x – Upgrades)

SHEET 18A OF 18


Contact Function Allocation			
Pin	Twisted with	Description	Comments
E1	n/a	Open	Not used
E2	n/a	Open	Not used
E3	n/a	Open	Not used
E4	n/a	Open	Not used
E5	n/a	Open	Not used
E6	n/a	Open	Not used
F1	G1	MC bias 1	MC T/R Bias Line 17 for "A", 25 for "B"
F2	G2	MC bias 2	MC T/R Bias Line 18 for "A", 26 for "B"
F3	G3	MC bias 3	MC T/R Bias Line 19 for "A", 27 for "B"
F4	G4	MC bias 4	MC T/R Bias Line 20 for "A", 28 for "B"
F5	G5	MC bias 5	MC T/R Bias Line 21 for "A", 29 for "B"
F6	G6	MC bias 6	MC T/R Bias Line 22 for "A", 30 for "B"
F7	G7	MC bias 7	MC T/R Bias Line 23 for "A", 31 for "B"
F8	G8	MC bias 8	MC T/R Bias Line 24 for "A", 32 for "B"
F9	F10	body transmit enable	
F10	F9	body transmit enable return	
G1	F1	MC bias return 1	MC T/R Bias Line 17 for "A", 25 for "B"
G2	F2	MC bias return 2	MC T/R Bias Line 18 for "A", 26 for "B"
G3	F3	MC bias return 3	MC T/R Bias Line 19 for "A", 27 for "B"
G4	F4	MC bias return 4	MC T/R Bias Line 20 for "A", 28 for "B"
G5	F5	MC bias return 5	MC T/R Bias Line 21 for "A", 29 for "B"
G6	F6	MC bias return 6	MC T/R Bias Line 22 for "A", 30 for "B"
G7	F7	MC bias return 7	MC T/R Bias Line 23 for "A", 31 for "B"
G8	F8	MC bias return 8	MC T/R Bias Line 24 for "A", 32 for "B"
G9	n/a	Reserved	
G10	n/a	Reserved	
H1	n/a	Reserved	
H2	n/a	Reserved	
H3	n/a	Reserved	
H4	n/a	Reserved	
H5	n/a	Reserved	
H6	n/a	Reserved	
H7	n/a	Reserved	
H8	n/a	Reserved	
H9	n/a	Reserved	
H10	n/a	Reserved	
I1	n/a	Reserved	
I2	I5	+15Vdc (auxiliary & preamp power)	
I3	I6	+10Vdc (auxiliary & preamp power)	
I4	n/a	-15Vdc (auxiliary power)	
I5	I2	Return for aux power	
I6	I3	Return for aux power	
I7	I8	coil ID (Dallas OneWire protocol)	All coils MUST have coil ID.
I8	I7	coil ID return (Dallas OneWire protocol)	All coils MUST have coil ID.
I9	I10	coil present'	
I10	I9	coil present return'	

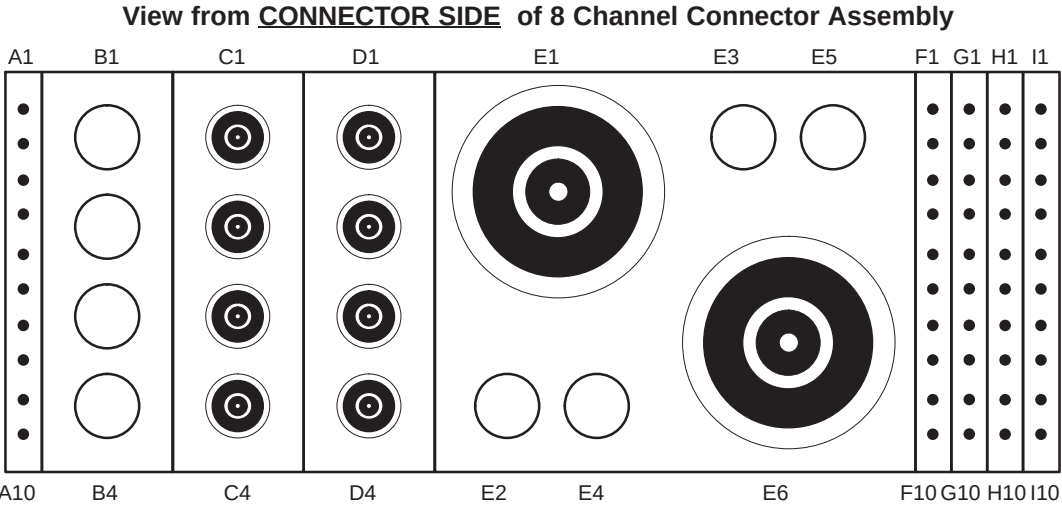

Contact Function Allocation			
Pin	Twisted with	Description	Comments
A1	A2	coil present	
A2	A1	coil present return	
A3	A4	multiplexerSelect0	RS-422+
A4	A3	multiplexerSelect0 return	RS-422-
A5	A6	multiplexerSelect1	RS-422+
A6	A5	multiplexerSelect1 return	RS-422-
A7	A8	multiplexerSelect2	RS-422+
A8	A7	multiplexerSelect2 return	RS-422-
A9	A10	multiplexerSelect3	RS-422+
A10	A9	multiplexerSelect3 return	RS-422-
B1	n/a		Not used
B2	n/a		Not used
B3	n/a		Not used
B4	n/a		Not used
C1	n/a	MC-1A,B or MC-9B	Preamp out 1, no bias, DC blocked
C2	n/a	MC-2A,B or MC-10B	Preamp out 2, no bias, DC blocked
C3	n/a	MC-3A,B or MC-11B	Preamp out 3, no bias, DC blocked
C4	n/a	MC-4A,B or MC-12B	Preamp out 4, no bias, DC blocked
D1	n/a	MC-5A,B or MC-13B	Preamp out 5, no bias, DC blocked
D2	n/a	MC-6A,B or MC-14B	Preamp out 6, no bias, DC blocked
D3	n/a	MC-7A,B or MC-15B	Preamp out 7, no bias, DC blocked
D4	n/a	MC-8A,B or MC-16B	Preamp out 8, no bias, DC blocked

G3 HDx (14.x) & HDxt (15.x) MGD/RRF/RF SUBSYSTEM
MCR III Tool – 8 Channel connector
(14.x – Forward Production)
SHEET 18B OF 18

Contact Function Allocation			
Pin	Twisted with	Description	Comments
E1	n/a	Open	Not used
E2	n/a	Open	Not used
E3	n/a	Open	Not used
E4	n/a	Open	Not used
E5	n/a	Open	Not used
E6	n/a	Open	Not used
F1	G1	MC bias 1	MC T/R Bias Line 17 for "A", 25 for "B"
F2	G2	MC bias 2	MC T/R Bias Line 18 for "A", 26 for "B"
F3	G3	MC bias 3	MC T/R Bias Line 19 for "A", 27 for "B"
F4	G4	MC bias 4	MC T/R Bias Line 20 for "A", 28 for "B"
F5	G5	MC bias 5	MC T/R Bias Line 21 for "A", 29 for "B"
F6	G6	MC bias 6	MC T/R Bias Line 22 for "A", 30 for "B"
F7	G7	MC bias 7	MC T/R Bias Line 23 for "A", 31 for "B"
F8	G8	MC bias 8	MC T/R Bias Line 24 for "A", 32 for "B"
F9	F10	body transmit enable	
F10	F9	body transmit enable return	
G1	F1	MC bias return 1	MC T/R Bias Line 17 for "A", 25 for "B"
G2	F2	MC bias return 2	MC T/R Bias Line 18 for "A", 26 for "B"
G3	F3	MC bias return 3	MC T/R Bias Line 19 for "A", 27 for "B"
G4	F4	MC bias return 4	MC T/R Bias Line 20 for "A", 28 for "B"
G5	F5	MC bias return 5	MC T/R Bias Line 21 for "A", 29 for "B"
G6	F6	MC bias return 6	MC T/R Bias Line 22 for "A", 30 for "B"
G7	F7	MC bias return 7	MC T/R Bias Line 23 for "A", 31 for "B"
G8	F8	MC bias return 8	MC T/R Bias Line 24 for "A", 32 for "B"
G9	G10	Preamp Enable	
G10	G9	Preamp Enable Return	
H1	n/a	Reserved	
H2	n/a	Reserved	
H3	n/a	Reserved	
H4	n/a	Reserved	
H5	n/a	Reserved	
H6	n/a	Reserved	
H7	n/a	Reserved	
H8	n/a	Reserved	
H9	n/a	Reserved	
H10	n/a	Reserved	
I1	n/a	Reserved	
I2	I5	+15Vdc (auxiliary & preamp power)	
I3	I6	+10Vdc (auxiliary & preamp power)	
I4	n/a	-15Vdc (auxiliary power)	
I5	I2	Return for aux power	
I6	I3	Return for aux power	
I7	I8	coil ID (Dallas OneWire protocol)	All coils MUST have coil ID.
I8	I7	coil ID return (Dallas OneWire protocol)	All coils MUST have coil ID.
I9	I10	coil present'	
I10	I9	coil present return'	

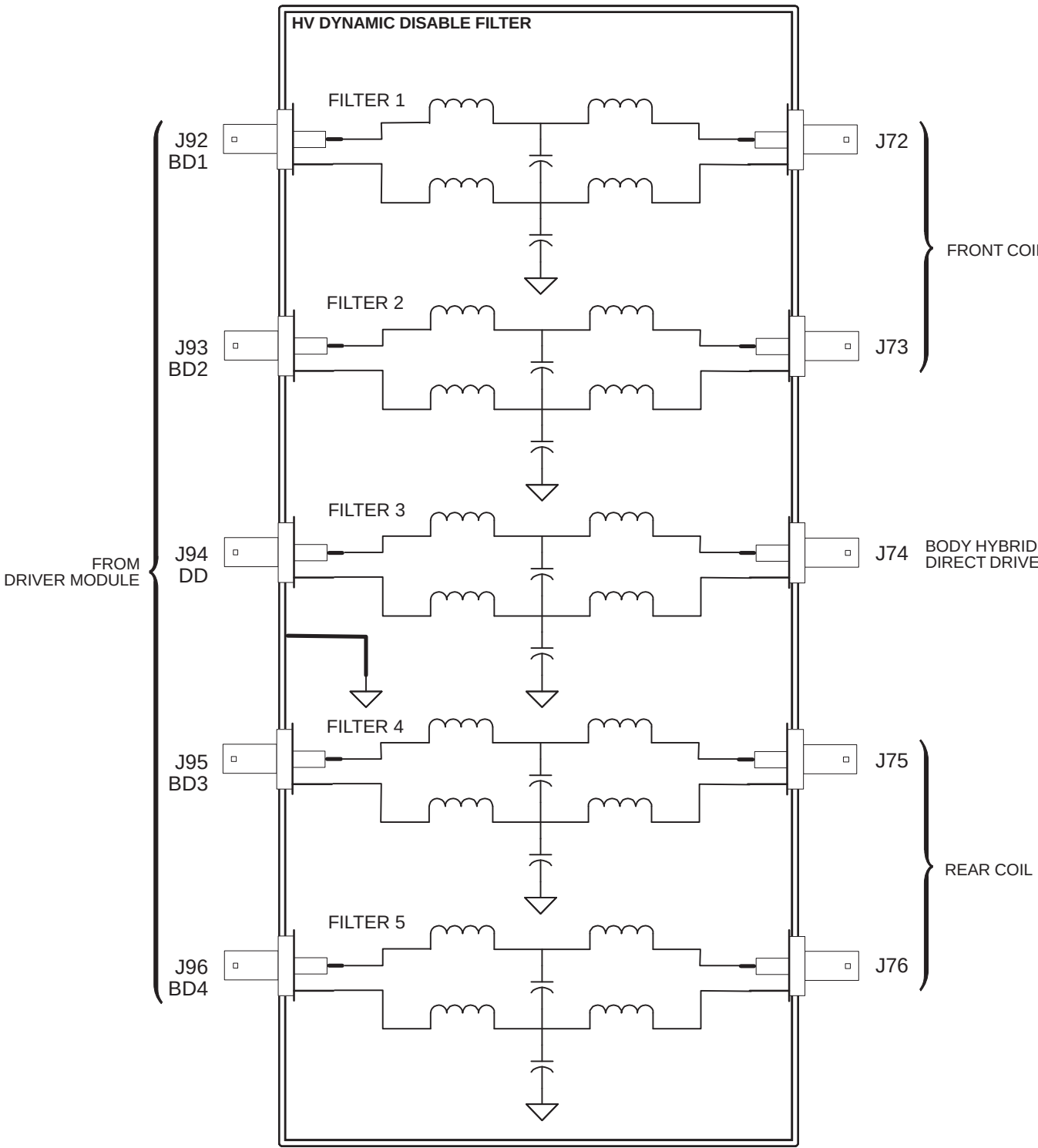

NOTE: Connector shields are electrically isolated from enclosure housing.

COOLING SYSTEM

TABLE OF CONTENTS

Water Chillers Subsystem Group Interconnect Diagrams

<u>ILLUSTRATION</u>	<u>PAGE 1-1</u>
1 OUTDOOR TSCC & RCP	
2 INDOOR WATER COOLED TSCC	
3 INDOOR AIR COOLED TSCC	
4 OUTDOOR TGWC & RCP	
5 INDOOR WATER COOLED TGWC	
6 INDOOR AIR COOLED TGWC	
7 DUAL INDOOR MRCC	
8 INDOOR MRCC	
9 DUAL OUTDOOR MRCC	
10 OUTDOOR MRCC	
11 INDOOR GWHX	

Note:

For all Chillers and Cooling System variations (SEE Vendor Manuals on MR Service Methods CDROM)

