

FCT SERVICE NOTE

SN05110812a


Date: November, 08, 2005

Subject: Tube Identification Certification Tool (TIC)

Effectivity: Irix/Pegasus systems with LS2002M5 software and **G.6 upgrade**
(FMI 25363)

Details:

1. With the release of Irix Patch G.6 via FMI 25363 (released FW44 – 2005) a new feature, SmarTube has been implemented. In addition to numerous quality and safety fixes, this software introduces SmarTube, which requires a Field Engineer to verify and identify of the tube installed on the system. If this selection is NOT made, OR, if the tube installed is a non-GE supplied tube, the following message will pop-up on each boot and at other times during operation. For more complete information, please see FMI 25363 doc. 5142331-100 or Tips & Workarounds 5115418-100.

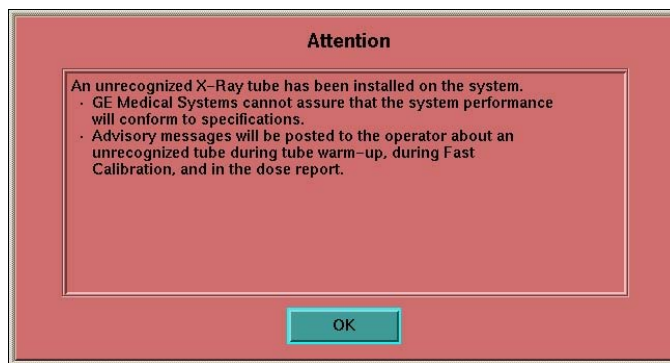

2. **The following service events will require the tube to be identified:**

- After x-ray tube is changed
- Upon Alignments
- After Load From Cold


Note: FE Service Key will be required to enable tube selection. If Service Key is not installed, a 30 day window will be provided for the FE to return to the site with the service key and make tube selection without adversely affecting the customer.

3. **Selection Of Tube Identity:** Follow steps listed below or Flowchart on next page!

Note: Follow all established safety procedures before entering gantry to verify tube!

1. Verify Tube identity by checking the tube insert model # on rating plate. For GE Performix Tube, model# is 2120785 or 2120785-2.
2. Open Common Service Desktop (CSD)
3. Install GE Class M (FE) service key.
4. Select Configuration
5. Select Tube Identification Certification (TIC) tool. (Do NOT double click on tool, that will open multiple windows with a possible error message “**service key information not available**”, if this happens close windows and retry, or reboot!
6. Make tube selection, (**‘GE Medical Systems Performix Tube’ or ‘All Others’**) in the TIC Tool. The tool is in Common Service Desk Top under Configuration.
7. Select “Confirm” then “Accept”!
8. Select “Dismiss” This will exit the tool.
9. Reboot Scanner.

Selection of Tube Identity

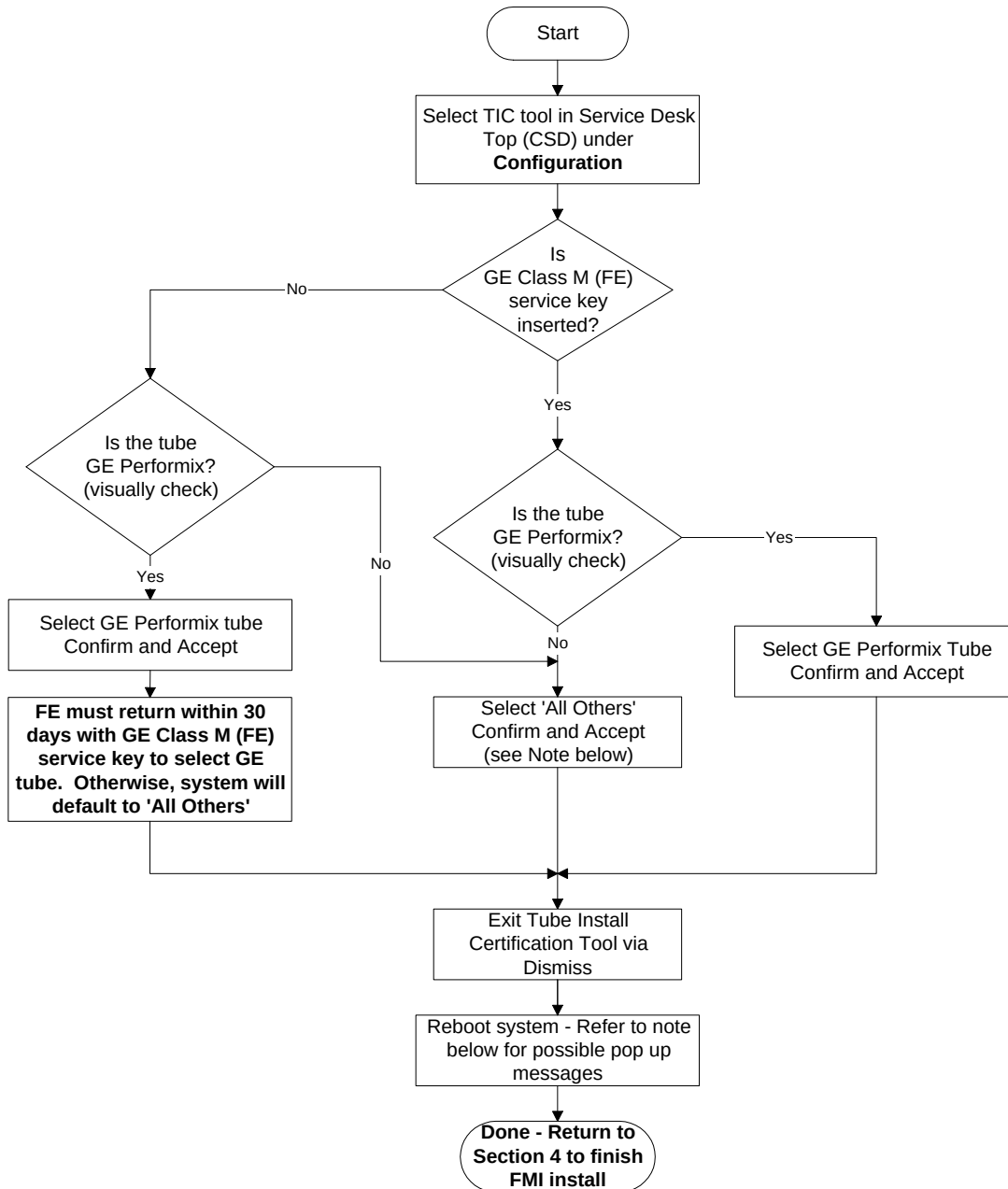

Note: Possible Pop Up messages at start up:

1. If the tube selection is GE tube without time limit, no pop up messages relating to tube identity on boot up.
2. If the tube selection is GE with time limit, a message with days remaining for GE FE to verify tube will be shown on boot up until GE FE removes the time limit with service key. If not verified by GE FE within the time limit, the system will default to 'All Others'.
3. If the tube is set to 'All Others', un-recognized x-ray tube installed message will always be shown at boot up.