

FCT SERVICE NOTE SN06101010A

Class C Service information enclosed.

Date: 10-Aug-06

Subject: IG, VIG, DARC2, VDARC Computer

Effectivity: LightSpeed 1x, 2x, 3x, 4x, 5x, 7x and HighSpeed QX/I GRE Host Console

Introduction:

The Intel Westville computer inside CT's IG, VIG, DARC2, and VDARC has gone end-of-life (EOL) and is replaced by the Intel Jarrell computer. GRE/Linux host console applications and diagnostics software have been updated via software release 06MW29.7 to support Jarrell IG, VIG, DARC2, and VDARC subassemblies.

Purpose:

The purpose of this service note is to introduce this EOL change to the field engineering team. It is to provide important information to support troubleshooting and optimize host console performance.

"Jarrell-based" subassemblies will begin to appear out in the field some time after October, 2006. Milwaukee and Hino manufacturing facilities are targeted to cut-in Jarrell subassemblies in October, 2006.

Prerequisites:

Software release 06MW29.7, or later, is required for "Jarrell-based" IG, VIG, DARC2, and VDARC to function. This software release supports Jarrell IG, VIG, DARC2, and VDARC functionality, applications, and diagnostics.

Software release 06MW29.7, software collector p/n 5181591, includes the following:

- 06MW29.7 Applications DVDROM,
- Linux OS 4.3.16 DVDROM, and
- Serial Over LAN for Jarrell and Westville-based DARC/VDARC CDROM.

How to Identify:

Westville and Jarrell subassemblies may be identified by their physical differences, part number differences, and via remote query using Linux command, "dmidecode."

Refer to the pictures below, illustrating their physical differences.

Refer to the **FRUs** section to for their part number labeling differences between Westville and "Jarrell."

Refer to the following steps for "dmidecode" instruction, how this Linux command may be used to remote query the host console to determine its configuration:

1. Login into the subassembly of interest (IG, VIG, DARC2, or VDARC) as "root"
2. Type **dmidecode | less**
3. View the "Version" line towards top of the response (refer to example response below)
 - o If **SWV25**, then WESTVILLE is installed
 - o If **SE7520JR23**, then JARRELL is installed

Westville Response (for IG, VIG, DARC2, VDARC):

NOTE: All use the same Intel Westville computer → response will be the same.

```
# dmidecode 2.6
# dmidecode 2.2
SMBIOS 2.3 present.
73 structures occupying 2782 bytes.
Table at 0x000F1260.
Handle 0x0000
DMI type 0, 20 bytes.
BIOS Information
  Vendor: Intel Corporation
  Version: SWV25.86B.0257.P42.0504011224
  Release Date: 04/01/2005
  ...etc....
```

Jarrell Response (for IG, VIG, DARC2, VDARC):

NOTE: All use the same Intel Jarrell computer → response will be the same.


```
# dmidecode 2.6
SMBIOS 2.3 present.
59 structures occupying 2536 bytes.
Table at 0x000FD4A0.
Handle 0x0000
DMI type 0, 20 bytes.
BIOS Information
  Vendor: Intel Corporation
  Version: SE7520JR23.86B.P.10.00.0087.120820051348
  Release Date: 12/08/2005
  ...etc....
```

Westville – Jarrell Physical Differences / Pictures

Westville DARC Picture -- Front




Jarrell DARC Picture -- Front

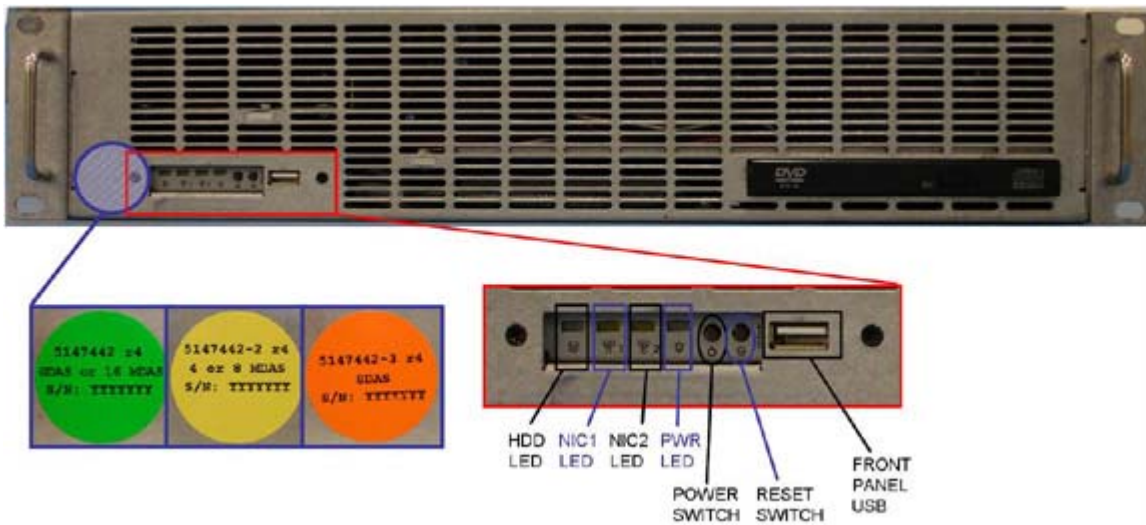

Westville IG -- Front



Jarrell IG -- Front

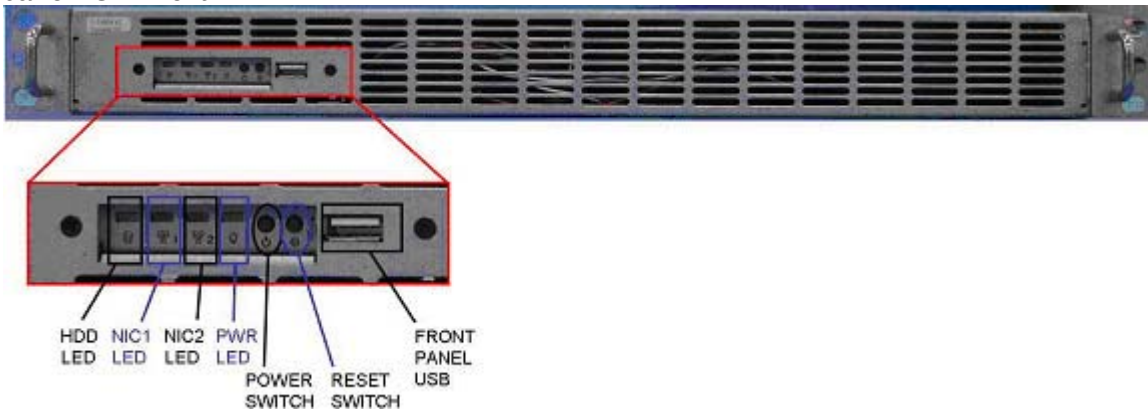

Training / Product Features:

Hardware

The Jarrell IG (p/n 5147443) requires a shielded Ethernet cable (p/n 5127222-2) to connect it to the DARC2 to maintain electromagnetic compatibility (EMC). "FE Assist Notes" have been added to FEMC to advise the FE of this requirement. This note is seen when the Jarrell IG is being ordered via FEMC. This shielded Ethernet cable already is present in the VCT / GOC5 configuration.

In multi-IG (or VIG) configurations, Westville and Jarrell IGs (or VIGs) may coexist -- mixed configurations are allowed. In a mixed configuration, the Jarrell IG (or VIG) must be in IG (VIG) Slot #1 for optimum image generation performance. IG (VIG) Slot #1 is the position closest to the top of the DARC2 (VDARC).

Software – Release 06MW29.7

Software release 06MW29.7, p/n 5181591, is required to support Jarrell-based IG, VIG, DARC2, and VDARC subassemblies. Refer to the **FRUs** section of this service note for these part numbers.

Host Console Diagnostics

GRE/host console diagnostics software has the same "look and feel," but has been expanded to test Jarrell-based IG, VIG, DARC2, and VDARC in software release 06MW29.7.

Service Toolkit:

For Class C/M customers / customer sites, software release 06MW29.7 requires the following Class C and M software:

- 06MW29.7 Class C p/n 5192235, distributed per TAB #182 distribution list
- 06MW29.7 Class M p/n 5181129, distributed per TAB #182 distribution list

Related Service Note(s)

SN06101011A, Introductory Service Note for Limited Access System Installations

SN06101012A, Introductory Service Note for VCT 64-Slice HALO DAS/Detector

Service Docs:

Service information CDROMs that ship in the technical publication collector with all systems have been updated to include Jarrell subassembly support. The following procedures / sections have been updated:

- Software Installation, 06MW29.7 Load from Cold
- Parts Replacement
- Illustrated Parts List
- GOCx Interconnect Diagram

FRUs:

The following field replaceable units (FRUs) have been created to support the Jarrell-based IG, VIG, DARC2, and VDARC:

P/N "Jarrell-based" FRU	Part Description	Qty / System	Replaces P/N
5147442	DARC 4/8/16slice GDAS & 16slice MDAS, Jarrell	1 1 1	5114772 2362870 2362870-10
5147442-2	DARC 4/8slice MDAS (console upgrade), Jarrell	1 1 1	5114772-2 2362870-2 2362870-12
5147442-3	DARC 4slice SDAS (console upgrade), Jarrell	1 1 1	5114772-3 2362870-3 2362870-13
5159833	VDARC w/VDIP, Jarrell	1	5114523-2
5147443	IG, Jarrell	1	2362872
5159834-2	VIG w/VRAC2, Jarrell	3	5114533-2

NOTE:

Software release 06MW29.7, p/n 5181591, or later, is required to support Jarrell-based IG, VIG, DARC2, and VDARC.

This software collector is orderable via the traditional CoakleyTech order site, <http://coakleytech.med.ge.com/cgi-win/default.aspx> for international field engineers. For USA field engineers, it is orderable via the GEMSbuy CoakleyTech Punchout site (GE Healthcare homepage → Healthcare Links → Tools → Self-Service Purchasing (SSP/GEMSbuy)).

Planned Maintenance:

No changes to planned maintenance schedules; no additions.

Mike Dejewski
LSI, MI & CT