

GE Healthcare Functional CT (FCT) Class M eService Notes are copyright (c) to General Electric Company Inc with all rights reserved. Class M eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare prohibited.
Duplication and distribution of any class of FCT eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

FCT SERVICE NOTE

SN06101806a

Class M information is Property of GE - For GE Service Personnel Only. Disclosure to Third Parties Prohibited.


Date: October 20, 2006

Subject: New Service Tool is now available: Introduction of the gesyslog File Viewer

Effectivity: CT and PET systems

Details: Current troubleshooting methodology involves gesyslog analysis. This can be a tedious and time consuming process of scanning through the error log (syslog) looking for a pattern of errors associated with the failure you are investigating.

Resolution: Introduction of the gesyslog File Viewer Service Tool. This tool can be loaded on a PC/Laptop or can be run directly from a CD media on any CT and PET system. The software tool is designed to search for patterns of key system failure events. When used, as a troubleshooting aid is designed to reduce the error log to a manageable "Set" of messages. It uses a color-coding scheme to provide and describe all the messages involved with key failure events. This reduces the time needed to determine which error is the root cause of the failure. Presently, This tool is designed for use on CT and PET systems only.

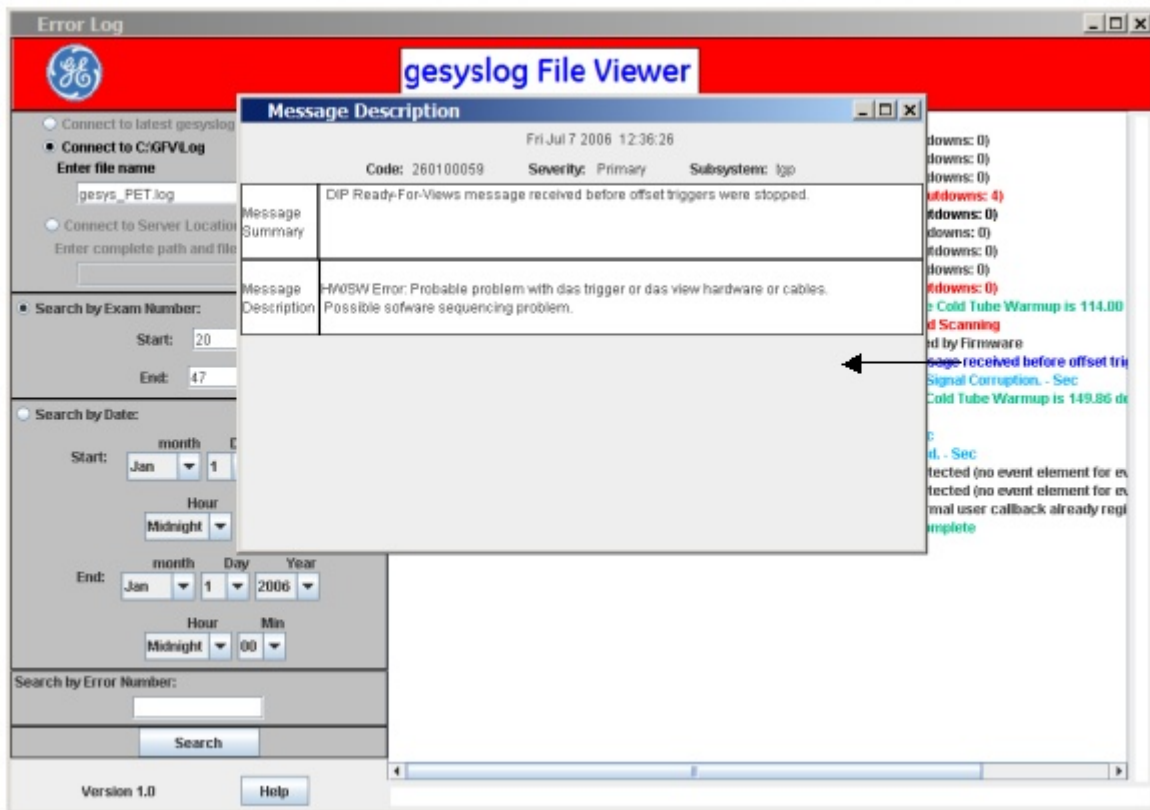

GE Healthcare Functional CT (FCT) Class M eService Notes are copyright (c) to General Electric Company Inc with all rights reserved. Class M eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare prohibited.
Duplication and distribution of any class of FCT eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

HOW TO GET THE gesyslog File Viewer Tool:

This tool will not be delivered or published on a CD to all Field Service personnel. It is ONLY available as a Download .ISO file at:

http://www.ct.med.ge.com/tech_pubs/fct/system/fct_home.htm

Scroll to the bottom of the Service Engineering Website:

1. Select: *Downloads*
2. Read the Proprietary warning statement, then Click on the blue lettering starting with:
"I certify that I am a GE Healthcare Employee...."
3. Select and download the file: **GFV_1.0.zip** to your computer.
4. Once downloaded unzip the file into a new (temp) folder.
5. Burn the .iso file to a blank CD media disk. Instructions on how to burn an .iso file is included with the zipped file.

Tutorial- Demo/Learn To Use: (Take 5 minutes to walk through the demo)

1. Load the software on your Service Laptop first, to become familiar with its use. If "Autorun" is enabled the software will guide you through the installation. If "autorun" is NOT enabled:
 - a. Select: Start
 - b. Select: Run
 - c. Browse to the media disk
 - d. Select: INSTALL.EXE
2. Open the software tool:
 - a. Select: Start
 - b. Select: Programs
 - c. Select: MI&CT gesyslog File Viewer folder
 - d. Select: MI&CT gesyslog File Viewer Program
3. From the main window, select the "Help" button. **Go to "Demo/Learn To Use"**.
4. Follow the Demonstration "Walk through" as an introduction tutorial. This should be sufficient to give you enough knowledge about the tools' operation to begin using it.

Updates will be announced periodically. You should check the Download Section on the MI&CT Service Website regularly for Service Documentation updates and Service Tool Downloads.

If you have problems running the software:

Locate the Readme.txt file on the CD media and read it for further instructions. Also, you can just select the "Help.htm" file found on the disk and your browser will open the Help.htm page. Refer to "Issues/Error messages" link.

Douglas R. Hofstetter
MI&CT Service Engineering