

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.
Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

MI&CT ESERVICE NOTE SN07072611a

Class C information is Property of GE.

This short-term informational note is invalid 6 months after the issue date.

Date: August 28, 2007

Subject: Rotating Index Lock Inspection prior to Gantry Installation

Effectivity: All Lightspeed Family (including VCT), Brightspeed and Brightspeed Select

Details:


Recently, we've received reports from the field that the rotational locking pin failed on the field during the installation or the first service. Root cause is the pins may got damaged or broken by the vibration and shock during the transportation. This leads the teeth could not be engaged fully. Corrective actions on manufacture procedure and supplier quality have been finished.

Resolution:

Reminder Notice: All FE's to follow the service manual strictly and verify the lock before any service actions.

To prevent un-commanded motion of the rotating assembly, the rotational locking pin must be engaged anytime the rotating assembly is serviced.

Rotational Lock Assembly

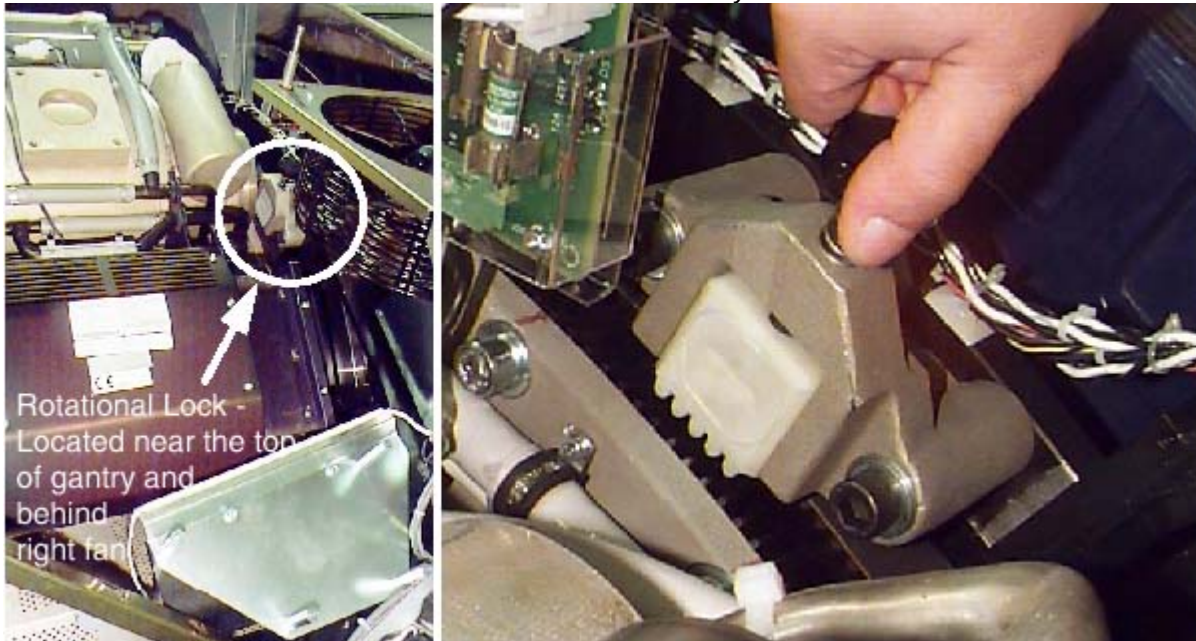

NOTE: For the 5.x Gantry, the fan shown in the illustration above is part of the top cover and is not on the gantry as shown.

The rotational lock is located on the rear side of the gantry, near the top. It is positioned directly across from the teeth in the rotational assembly.

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.


Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

To operate the lock:

Turn the handle clockwise until the teeth on the lock fully engage the teeth on the rotating assembly. If necessary, rock the rotating assembly slightly to align the teeth. Hand tighten until snug. Do not over-tighten. **Visually verify that the teeth are engaged and manually rotate the gantry to verify if it is locked successfully.**

Turn the handle counter-clockwise until the teeth on the lock and the rotating assembly are fully disengaged and the teeth clear each other sufficiently.

Rotational Lock Assembly Operation

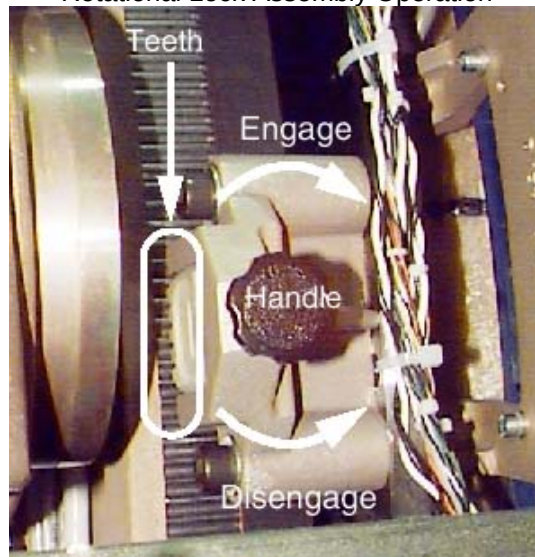

Du Ning
MI&CT Engineering, China

MI&CT Engineering Manager Reviewer: _____ Mark Nyman _____