

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

**MI&CT ESERVICE NOTE
SN07082410a**

Class C information is Property of GE.

This short-term informational note is invalid 6 months after the issue date.

Date: 08/24/2007

Subject: Global Table Stain on Side Leg Panel

Effectivity: VCT 64, VCT32 and RT16

Details:


It is commonly seen that the Flappers installed on Table Base Covers make blue stain on Table Side Leg Panels as Table moves up and down over a period of time.



Resolution:

Engineering decided to remove the Flapper from the Base Cover as it is installed on the cover just for cosmetic reason.

This change will be applied to the forward production as soon as supplier's preparation is completed. For installed base, remove the Flapper as needed. The Flapper can be removed by removing two screws.

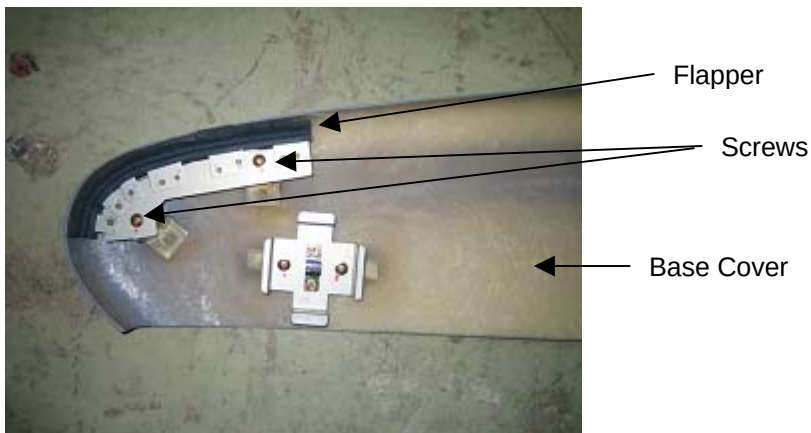

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.
Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

Affected Service Publication: N/A

PQR/SPR #: 13081401, 13116469

Takashi Uehara
Lead Service Integrator, GEYMS MI&CT Engineering

MI&CT Engineering Manager Signature_____Mark Nyman_____.