



***GE Medical Systems***

---

# **Technical Publication**

**Direction 2272430-100**

**Revision 0**

**Service Note**

**T-1432, Astec Power Supply  
Installation Instructions**

**RESTRICTED TO GENERAL ELECTRIC MEDICAL SYSTEMS (GEMS)**

THIS MATERIAL IS THE PROPERTY OF GEMS AND CONTAINS INFORMATION WHICH IS RESTRICTED TO USE BY GEMS MANUFACTURING, ENGINEERING AND SERVICE PERSONNEL. ANY REPRODUCTION IN WHOLE OR IN PART IS STRICTLY PROHIBITED. NEITHER THIS MATERIAL NOR ANY INFORMATION CONTAINED HEREIN SHALL BE DISCLOSED TO ANYONE OUTSIDE OF GEMS MANUFACTURING, ENGINEERING AND SERVICE PERSONNEL.



*Page Left Intentionally Blank*

## LEGAL NOTES

---

### GE EXCLUSIVE USE ONLY

Restricted To General Electric Medical Systems (GEMS)

This Material is the property of GEMS and contains information which is restricted to use by GEMS Manufacturing, Engineering and Service personnel. Any reproduction in whole or in part is strictly prohibited. Neither this material nor any information contained herein shall be disclosed to anyone outside of GEMS Manufacturing, Engineering and Service personnel.

If you are not a GE Medical Systems Employee, you must exit this material now.

### TRADEMARKS

All products and their name brands are trademarks of their respective holders.

### OMISSIONS & ERRORS

GE personnel, please use the GEMS CQA Process to report all omissions, errors, and defects in this publication.

### COPYRIGHTS

All Material Copyrighted(c) 2000 by the General Electric Company, All rights reserved.

### REVISION HISTORY

#### REVISION HISTORY

| Revision | Date     | Reason for change |
|----------|----------|-------------------|
| 0        | 11/06/00 | Initial Release   |
|          |          |                   |

# TABLE OF CONTENTS

|                                                                          |   |
|--------------------------------------------------------------------------|---|
| Legal Notes .....                                                        | 3 |
| GE EXCLUSIVE USE ONLY .....                                              | 3 |
| TRADEMARKS .....                                                         | 3 |
| OMISSIONS & ERRORS.....                                                  | 3 |
| COPYRIGHTS.....                                                          | 3 |
| REVISION HISTORY .....                                                   | 3 |
| Table of Contents.....                                                   | 4 |
| Power Supply Installation.....                                           | 5 |
| SECTION 1 – DIFFERENCES IN PACKAGING .....                               | 5 |
| SECTION 2 – INSTALLATION INSTRUCTIONS FOR NEW POWER SUPPLY (ASTEC) ..... | 5 |
| SECTION 3 – ADJUSTING THE POWER SUPPLY .....                             | 7 |

# POWER SUPPLY INSTALLATION


## SECTION 1 – DIFFERENCES IN PACKAGING

The new power supply by Astec is packaged with a sheet metal pan and the mounting screws necessary to attach it to the pan. The previous supply did not come with the pan or the hardware. This is necessary because the old pan will NOT support the new power supply. The new pan, however, will support both the new and old power supply.

## SECTION 2 – INSTALLATION INSTRUCTIONS FOR NEW POWER SUPPLY (ASTEC)

### Section 2.1 – Mounting the Power Supply to the Pan


When replacing the console power supply, remove the old power supply, including its pan. Then attach the new power supply to its pan using the screws shown in the drawing below.

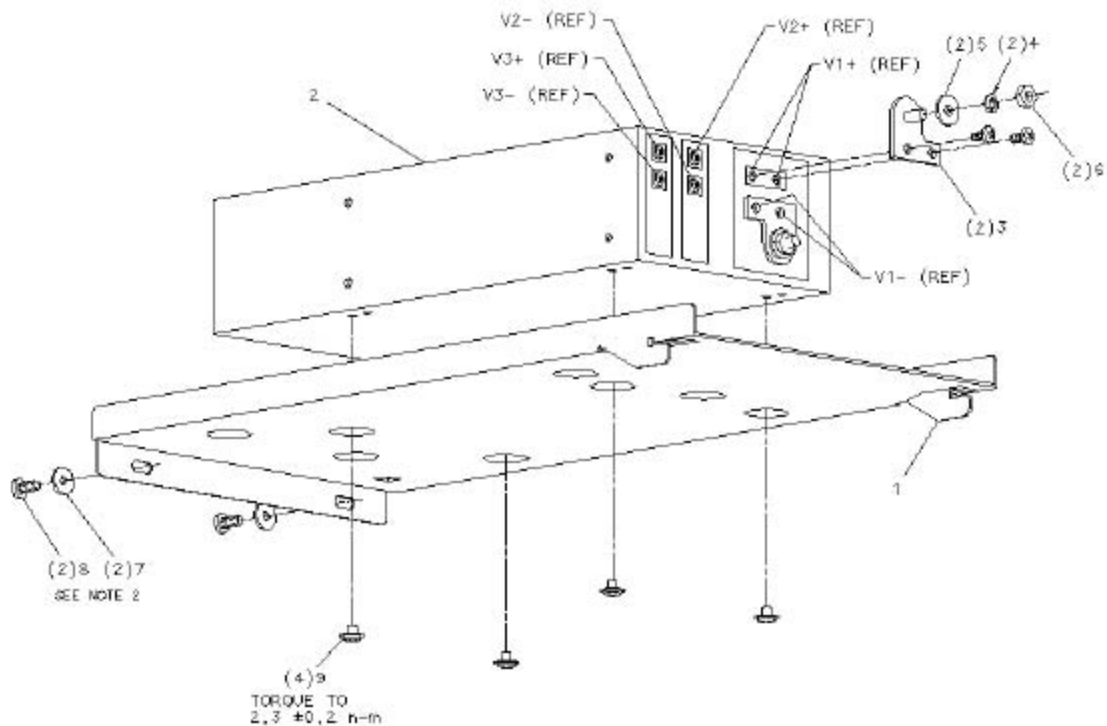

## Section 2.2 – Installing the Power Supply


When replacing an existing Power One power supply with the new Astec supply, use the diagram below to match wires to terminals. Each terminal in the new supply will have a corresponding wire, unlike the terminal strip in the Power One supply, which had several vacant positions.

### Power Supply DC Terminal Decoder { Power One Terminal numbers } to Astec terminal locations.

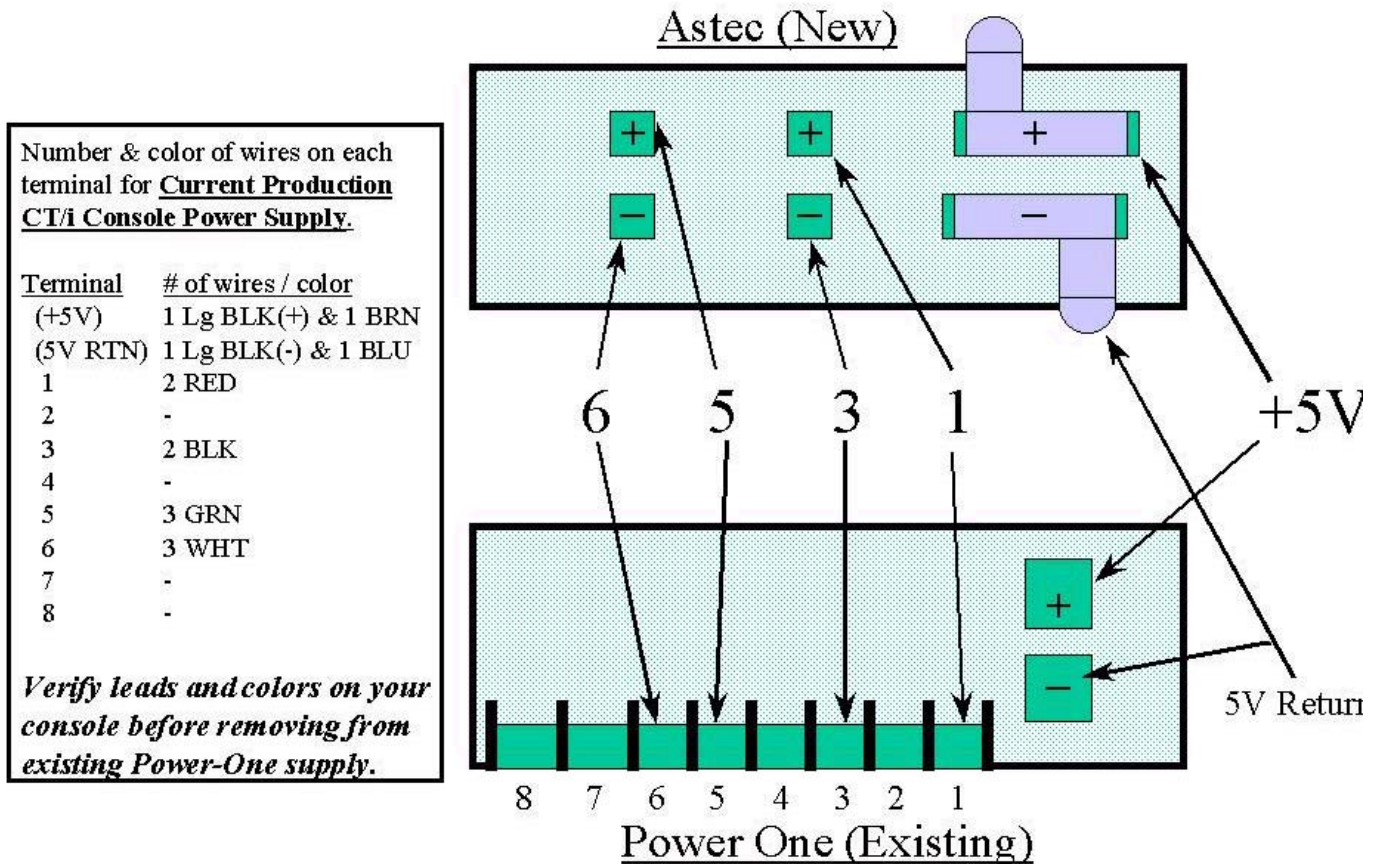

## SECTION 3 – ADJUSTING THE POWER SUPPLY

When replacing an existing Power One power supply with the new Astec supply, it will be necessary to adjust the Astec supply voltages. With the power supply under load and monitoring the power supply output, adjust the associated potentiometer for the readings listed in following diagram.

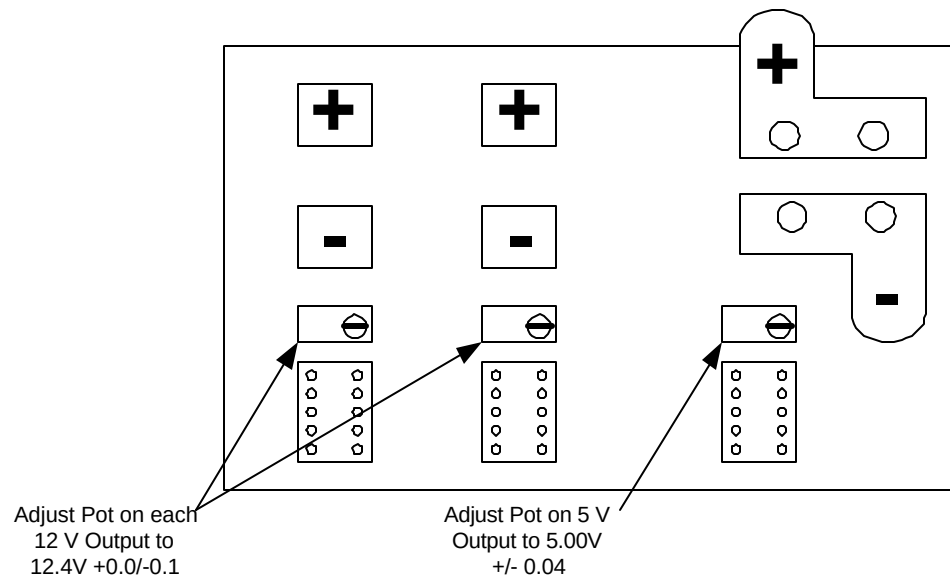

**GE Exclusive Use Only**

Direction 2272430-100, Rev 0

Service Note

**T-1432**, Astec Power Supply  
Installation Instrucionts

---



---

***GE Medical  
Systems***

*P.O. BOX 414*

*MILWAUKEE WISCONSIN  
53201-0414 U.S.A.*