


820-0241 – Install Instructions for Cryo Filter Kit (P/N: 205-0121-13) GE P/N: 5269685

Tools Required:

- 7/8" wrench or slip wrench
- Loctite[™] 565

Removal Procedure:

1. Perform shutdown of system in accordance with Lytron manual #820-0190.
2. Lockout/Tag out system in accordance with GE procedures.
3. Drain system. Drain system in accordance with Lytron manual #820-0190. Alternate draining through the pump housing is permissible.
4. Close ball valves on both sides of the filter.

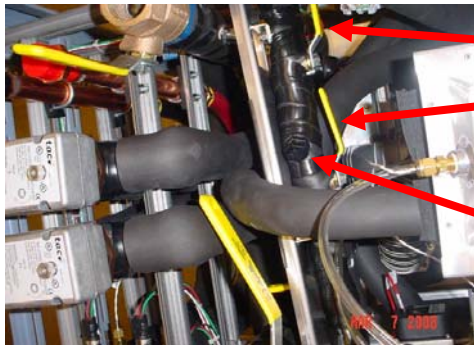

Isolation Ball Valves

Filter Retaining Nut

Photo 1

5. Using a 7/8" wrench or a Slip wrench, remove the filter retaining nut.
6. The filter is held to the bottom of the nut by friction – be careful when removing.
(Refer to Photo 2)

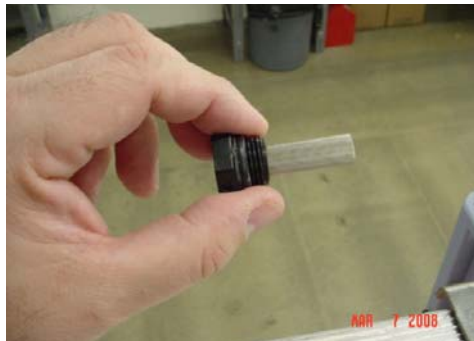


Photo 2

7. Clean or replace the filter screen.
8. Apply adequate amount of Loctite[™] 565 to the threaded portion of the retaining nut.
9. Place clean or new screen into the bottom of the retaining nut ensuring it is well seated.
10. Re-insert by holding the retaining nut and inserting entire assembly into manifold.
11. Tighten nut with 7/8" wrench or Slip wrench.

12. Re-open the two isolation valves you closed earlier.
13. Refill the system in accordance with Lytron manual #820-0190.
14. Start system and check for leaks.

End of Procedure