GE Healthcare

Discovery MR450, MR750 &
Optima MR450w
Pre-Magnet Delivery System Installation
Technical Publication

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File a report with

- Call 1–800–548–3366 and use option 8.
- Contact your local service coordinator for more information on this process.

Rev. 08/15/2003

© 2008–2009 General Electric Company
GE Medical Systems, a General Electric Company, going to market as GE Healthcare
3000 N. Grandview Boulevard
Waukesha, Wisconsin 53188
USA

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WARNING

(EN)

- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR, OR PATIENT FROM ELECTRIC SHOCK, OR FROM MECHANICAL OR OTHER HAZARDS.

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This page left blank intentionally.

REVISION HISTORY

REV	DATE	PRIMARY REASON FOR CHANGE
-----	------	---------------------------

1	May 12, 2008	Preliminary release for verification.
2	May 15, 2008	Section 1: Minor changes to wording. Section 2: Corrected spelling errors. Added power cord service loop in HEC cabinet. Release for validation.
3	Apr 17, 2009	Updated document title to Discovery MR450 & MR750 Pre-Magnet Delivery System Installation. Section 2-6: Added illustrations showing service loop on top of power box and the power wire connections to the terminal strip inside the power box (PQR's 13233956, 13236584). Section 2-7-1: Revised section to clarify the orientation of the SPW and Penetration Panel frames. Section 2-9: Added instructions to refer to Direction 5500102 for instructions on installing the curtain covers in the magnet room (PQR 13224769).
4	Nov 2, 2009	Updated document title to Discovery MR450, MR750 & MR450w Pre-Magnet Delivery System Installation. Section 2-7-1: Updated instructions for installing PEN and SPW mounting frames (PQR 13280034).

LIST OF EFFECTIVE PAGES

PAGE	REV	PAGE	REV	PAGE	REV	PAGE	REV	PAGE	REV
Title Page	4	TAB 2							
Damage in Trans.	–	COMPONENT							
Language Policy	4	INSTALLATION							
A	4	2-1 to 2-18	4						
i	4								
TAB 1		TAB 3							
GETTING STARTED		COMPLETION							
1-1 to 1-10	4	CHECKLIST							
		3-1 to 3-2	4						

TABLE OF CONTENTS

<u>SECTION</u>	<u>PAGE</u>
SECTION 1 – GETTING STARTED	1-1
SECTION 2 – COMPONENT INSTALLATION	2-1
SECTION 3 – COMPLETION CHECKLIST	3-1

SECTION 1 – GETTING STARTED

TABLE OF CONTENTS

<u>SECTION</u>	<u>PAGE</u>
1-1 INTRODUCTION	1-2
1-2 SITE READY CHECK FOR MECHANICAL INSTALLATION	1-2
1-3 NON-ENGLISH LABEL INSTALLATION	1-2
1-4 REQUIRED TOOLS	1-3
1-5 BASIC SYSTEM	1-4
1-6 PRODUCT DELIVERY INSTRUCTIONS	1-6
1-7 DAMAGE IN TRANSPORTATION	1-6
1-8 PRODUCT LOCATOR	1-7
1-9 INSTALLATION PROCEDURE	1-8
1-9-1 INSTALLATION FLOW	1-8
1-9-2 INSTALLATION FLOWCHART EXPLANATION	1-8
1-10 INSTALLATION FLOWCHART	1-9
1-10-1 PREREQUISITES FOR SYSTEM INSTALLATION	1-9
1-10-2 PRE-MAGNET DELIVERY MECHANICAL INSTALLATION	1-10

1-1 INTRODUCTION

This manual provides instructions to install the following components of the MR750 3.0T, MR450 1.5T or the MR450w System prior to delivery and installation of magnet.

- Main Disconnect Panel
- HEC Cabinet containing the Cooling Loops, and Gradient Blower.
- SPW (Second Penetration Wall) containing Gradient Filter, Scan Room Filters & I/F, and Ramp/Shim access
- HEC to Cryogen Chiller hoses and insulation
- SPW and PEN Panel mounting frames
- Magnet room Curtain Rail

This Direction covers the following catalog:

M7000BM – Preinstallation Collector for MR750 3.0T and MR450 1.5T

1-2 SITE READY CHECK FOR MECHANICAL INSTALLATION

Before equipment is delivered, one of the following applicable manuals must be complete to avoid delays and confusion.

- ☐ Pre-installation work must be complete. Refer to Direction 5500101, *Discovery MR750 3.0T Preinstallation manual*
- ☐ Pre-installation work must be complete. Refer to Direction 5500109, *Discovery MR450 1.5T Preinstallation manual*
- ☐ Pre-installation work must be complete. Refer to Direction 5670001, *Optima MR450w 1.5T Preinstallation manual*

1-3 NON-ENGLISH LABEL INSTALLATION

Components and accessories delivered to this site may have English labels attached to them. For non-English sites, optional multi-language labels have also been supplied.

If required, non-English labels should be attached over the English labels.

- For Surface Coils and Coil Cables – Apply the appropriate label over the English label. The label is included with the coil component.
- For all other labels, follow the label update procedures in the manual or on the applicable multi-language label sheet.

Verify at the end of the installation that the original English labels are covered by new labels representing the native language of the customer.

1-4 REQUIRED TOOLS**WARNING!**

FERROUS MATERIAL HAZARD! IF MAGNET IS ENERGIZED, TOOLS AND PARTS REQUIRED FOR THIS INSTALLATION THAT CONTAIN FERROUS MATERIAL AND WILL BE STRONGLY ATTRACTED TO MAGNET AND MAY BECOME DANGEROUS PROJECTILES. KEEP ALL FERROUS TOOLS AT LEAST 10 FEET AWAY FROM THE MAGNET.

Tools to install the system are listed below in Table 1-1.

TABLE 1-1
INSTALLATION EQUIPMENT


ITEM	GE PART NUMBER	DESCRIPTION
1	46-271138G1	Restricted Access Control Kit. Contains two plastic warning signs for posting at site during installation and service activity.
2	2319156	Aluminum platform ladder, 47.5 inches (1206.5mm) (See Note 1)
3	46-320273G3 or G4	Non-Magnetic Tool Kit – Universal (See Note 2) Both metric and inch Non-Magnetic Tool Kits needed. May substitute both of the following kits: • 46-320273G1 Non-Magnetic Tool Kit – Metric • 46-320273G2 Non-Magnetic Tool Kit – Inch
4	—	Volt Meter
5	—	Extension cords, with ground conductor
6	—	Power strip, grounded type, with minimum of five outlets
7	—	Plastic or aluminum flashlight
8	—	Non-magnetic level
9	—	Non-magnetic tape rule, 12 ft
10	—	Assorted drill bits
11	—	Inspection mirror
12	—	Hobby and utility knives
Note		
1	Supplied as part of Signa.	
2	Supplied by GE until turnover of system to customer, then available as part of a GE Cryogen and/or Service Contract.	

1-5 BASIC SYSTEM

The complete basic system, Illustration 1-1, consists of the following major equipment:

- Magnet Enclosure and Magnet Accessories
 - 1.5T (15 kilogauss) LCC Magnet
or
 - 3.0T (30 kilogauss) LCC300 Magnet
- Shield/Cryo Cooler Compressor Cabinet
- Gradient (XRMB) and RF body coils and Head Coil
- Main Disconnect Panel
- System electronics cabinets:
 - PGR Cabinet containing HFD Gradient drivers and Power Distribution Unit module with unregulated transformer 200/380/400/415/480 Volt, 50/60 Hz with power filter, CAM Chassis with boards for either 16 or 32 channels, components for the Volume Reconstruction Engine (VRE), Narrow Band RF Amplifier
 - HEC Cabinet containing the Cryo Coolers, Cooling Loops, and Gradient Blower.
 - PEN Cabinet containing the PHPS, Scan Room power supply, CAN Fiber I/F, T/R DD Driver module, and HO Shim Filter (optional)
 - SPW (Second Penetration Wall) containing Gradient Filter, Scan Room Filters & I/F, and Ramp/Shim access
- Operator Workspace equipment: GOC Computer Cabinet (containing Linux PC and GOCAA), LCD Widescreen monitor, keyboard, SCIM, SCSI Tower, mouse and mouse pad, and chair
- Pneumatic Patient Alert System
- Patient Transport Table and cradle
- Patient accessories such as: a phantom kit, patient log book, head cushion and sponges, chin and forehead straps, body wedges, knee cushions, and security/restraint straps
- Gating accessories which include: patient cardiac leads, peripheral gating probe, and respiratory bellows

1-5 BASIC SYSTEM (Continued)

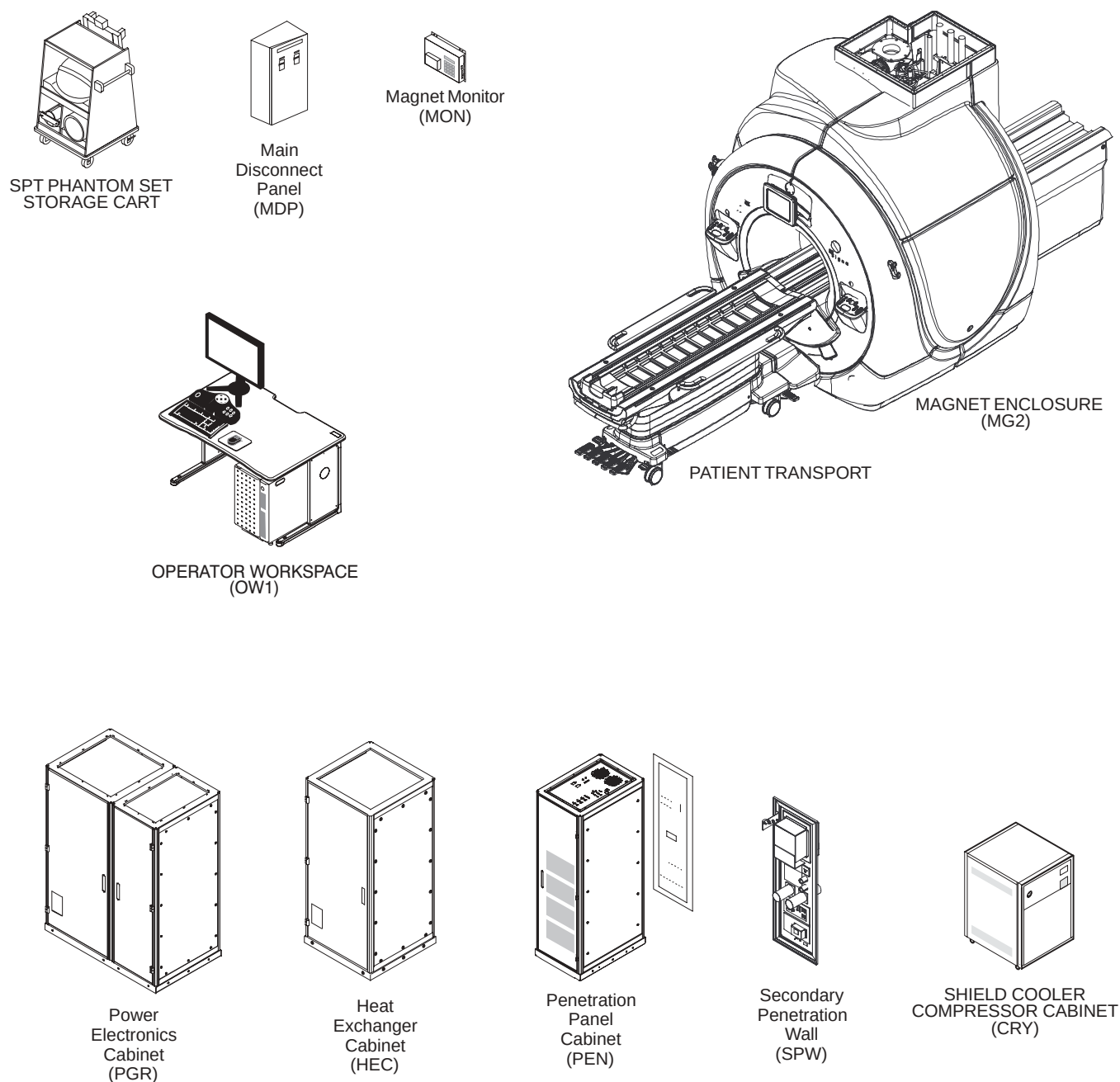SIGNA HDx 3.0T SYSTEM
ILLUSTRATION 1-1

1-6 PRODUCT DELIVERY INSTRUCTIONS

Note

Two or more people are required to accept delivery of equipment due to weight and size of cabinets and other components.

The “Shipping Document” lists catalog numbers delivered. Review and confirm that order is delivered complete. Check impact on installation schedule if Catalogs and/or packing boxes are missing and/or noted as shipped short.

Labels that are attached to the outside of packing boxes summarize box contents. The labels are color coded for different installation areas as follows:

- Green – Equipment Room
- Purple – Operator’s Room
- Orange – Magnet Room
- Yellow – Post-installation (Phantoms, Coils).

“Product Delivery Instructions” (PDI) specify box contents, part numbers, and shipping procedure. The PDI is numbered according to the catalog number (for example, PDI – M3333TA is for the Fixed Site Installation Kit). Lists of items included with each box are detailed by separate checklists, or a separate sheet that provides a summary of that box contents. Refer to PDI and packing lists for information specific to your shipment.

A set of service and operator manuals is delivered with the Fixed Site Kit. Refer to checklists packed with “Technical Publication” boxes for a list of delivered documents.

1-7 DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation “**damage in shipment**” written on **all** copies of the freight or express bill **before** delivery is accepted or “signed for” by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

File a report with

- Call 1-800-548-3366 and use option 8.
- Contact your local service coordinator for more information on this process.

1-8 PRODUCT LOCATOR

The Global Install Base Database (also known as Product Locator System) tracks shipment, trans-shipment, and field location of the serialized models. There are now two methods for submitting Product Locator information.

1. At this time, for **U.S. ONLY**, the preferred method of submitting information is the *FE Site Verification Web Site* at:
<http://gein2.med.ge.com/gib>

The FE Site Verification consists of three components that are available on the web from the main menu. They are:

- Install/deinstall product locator model and serial numbers
- Add/modify ship to address information
- Update CARES FE data for primary/secondary FEs

To obtain a copy of the FE training tutorial for using this website, a downloadable copy is available at the following:
Product Locator Support Central Page – <http://supportcentral.ge.com/15563>

2. One “Shipping Card” is filled out and submitted when shipped (extra cards are supplied for trans-shipments between storage and distribution points), and the “Installation Card” and extra shipment cards are attached.

Verify that serial and model number on each rating plate matches installation card numbers before removing installation card. Note that there may be one or more shipment cards and bar code labels with the installation card. These shipment cards are used to trace the transfer of serialized units between various inventory storage and distribution points until the product reaches its final installation destination. Process just the installation card and discard any extra shipment cards and labels.

1-9 INSTALLATION PROCEDURE

1-9-1 Installation Flow

Note

All on-site construction must be completed before equipment is delivered and installation starts. Attempting to install the system while construction is being completed will impact installation efficiency and further delay site completion. Making sure that all preinstallation and construction work is completed before equipment is delivered will usually result in an earlier turnover date.

The flow chart on the next page should be followed for an orderly and efficient installation of the Pre-magnet delivery components of the system being installed. Note that many procedures may be performed in parallel and may be performed in any order according to the specific situation of each site.

Note


The flow chart references the Tab number and/or Section number for the respective installation procedures. These procedures are similar to cable maps in that they are to be removed from this binder and used at the location where hardware installation is being performed. The Procedures must be returned to the binder for future reference upon completion of the installation.

This installation flowchart has been developed assuming that all components for this part of the installation have been delivered together. Make sure that every part required is available before starting.


1-9-2 Installation Flowchart Explanation

Shown below are three examples of the flow chart Installation Procedures and an explanation of the areas within the procedures.

Use the correct Magnet Manual that applies to the System being installed to locate the installation procedure listed below

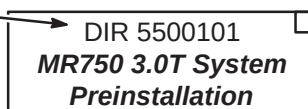
Tab within this manual where the Section listed below is located

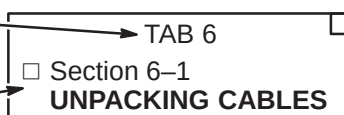Section within the Tab where the Installation Procedure is located



LEGEND FOR FOLLOWING INSTALLATION FLOWCHART

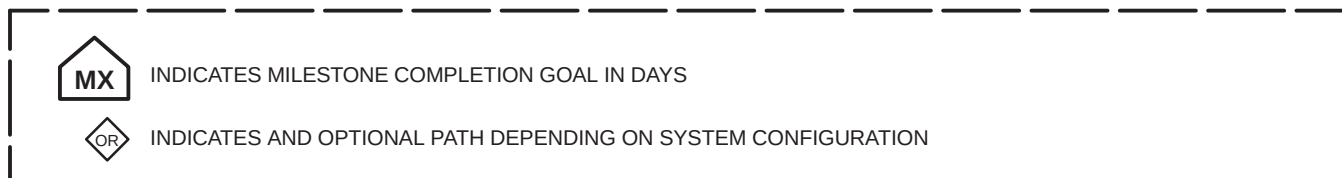
1-10 INSTALLATION FLOWCHART

1-10-1 PREREQUISITES FOR SYSTEM INSTALLATION

Site Installations Required Prior to System Mechanical Installation

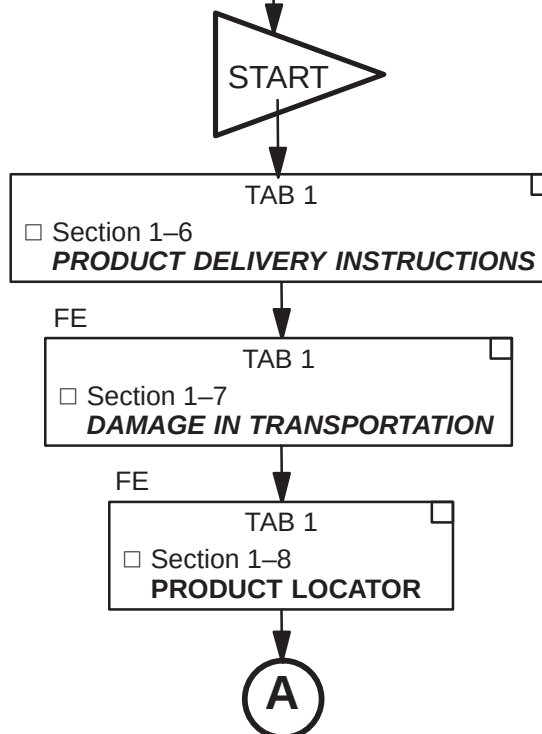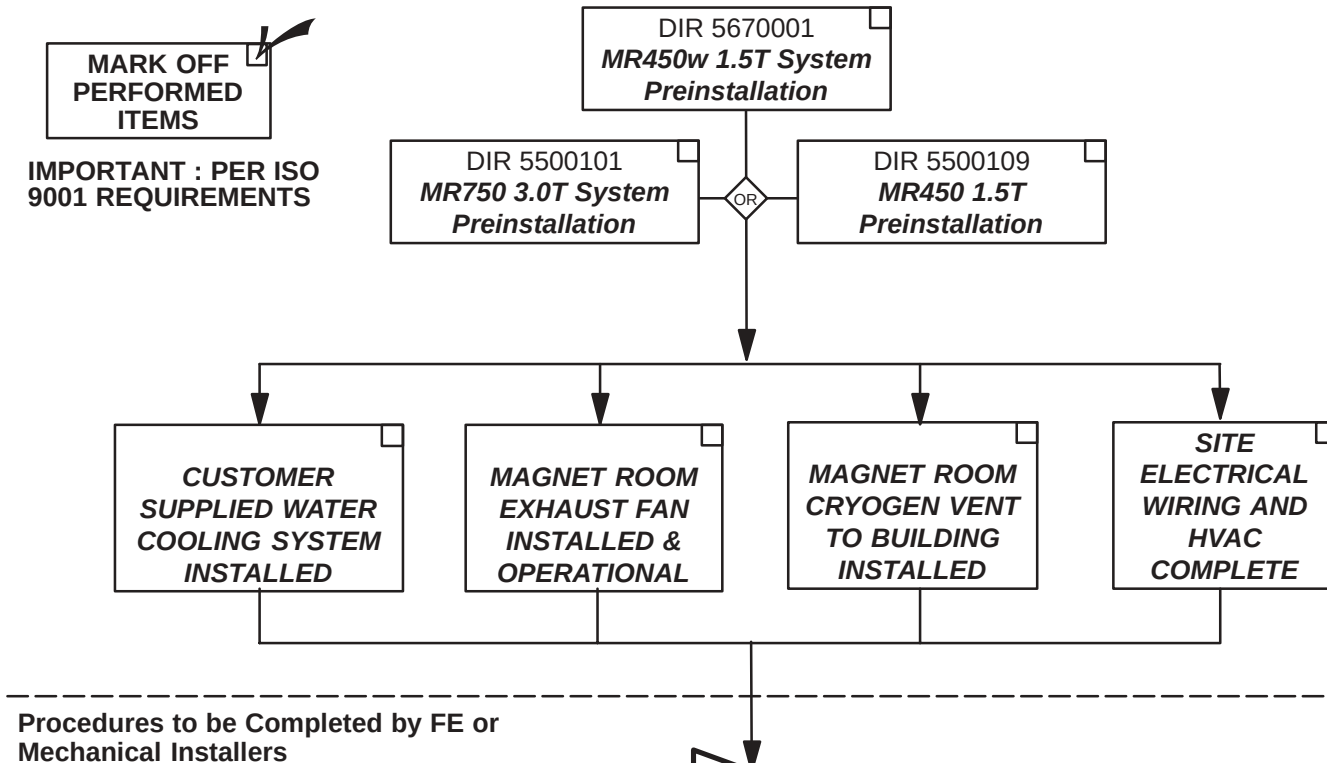1-10-2 PRE-MAGNET DELIVERY MECHANICAL INSTALLATION

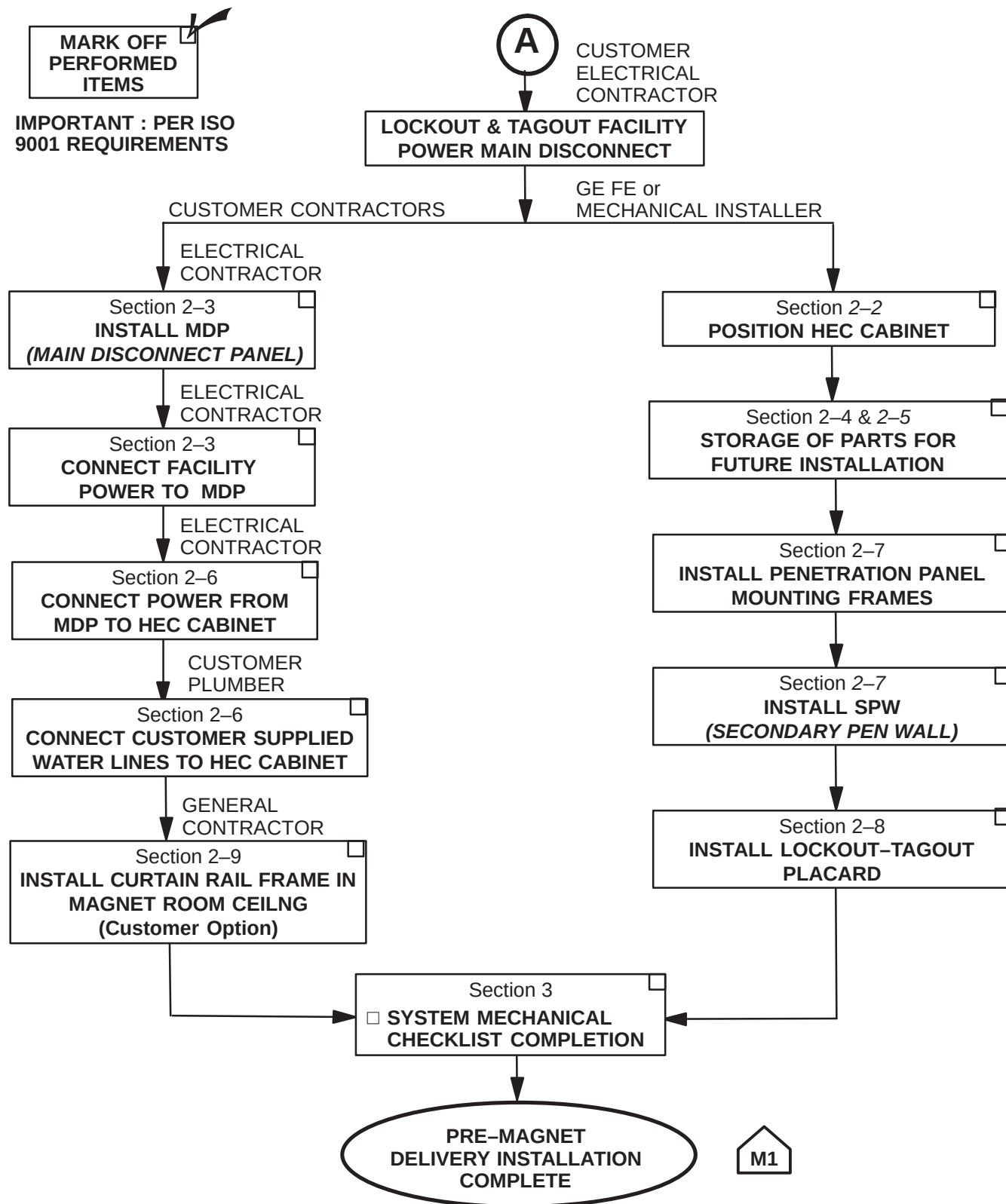SECTION 2 – COMPONENT INSTALLATION

TABLE OF CONTENTS

<u>SECTION</u>	<u>PAGE</u>
2-1 COMPONENT INSTALLATION REQUIREMENTS	2-2
2-2 HEC CABINET INSTALLATION	2-3
2-3 INSTALLATION OF MAIN DISCONNECT PANEL	2-8
2-4 HEC TO CRYO-COOLER HOSE KIT	2-8
2-5 HEC TO CRYO-COOLER HOSE INSULATION KIT	2-8
2-6 CUSTOMER CONNECTIONS TO HEC CABINET	2-9
2-7 INSTALL RF SCREEN ROOM WALL MOUNTING FRAMES AND SPW	2-10
2-8 INSTALL LOCK-OUT TAG-OUT PLACARD	2-12
2-9 INSTALL MAGNET ENCLOSURE CABLE CURTAIN (Customer Option for MR750 only) ...	2-13

2-1 COMPONENT INSTALLATION REQUIREMENTS

Positioning of cabinets and components is critical for installing cables and post-installation system operation. Refer to architectural site layouts for:

- Location of all cabinets and related equipment.
- For seismic anchoring (where required)

Equipment room cabinets must not be installed in an area with magnetic fields higher than 50 gauss. Refer to site layout plans for positioning of cabinets.

DUE TO WEIGHT OF CABINETS, AT LEAST TWO PEOPLE ARE REQUIRED TO MOVE ANY CABINET

2-2 HEC CABINET INSTALLATION

2-2-1 HEC CABINET POSITIONING PREPARATION

The illustration below shows a typical layout for the equipment room.

Refer to architectural site layout drawings to determine location of HEC cabinet.

- If available, use the cabinet floor layout templates for positioning cabinets
- If templates are not available, using the minimum spacing distances as noted below, apply tape on floor to mark final location of front corners of cabinets.


The PEN and PGR cabinets will be delivered and installed after installation of the magnet.

Legend:

A = 4 inch (102mm) minimum space distance from wall


B = 15 inch (380mm) minimum space distance at side of cabinet

 = Final location of PEN and PGR Cabinets. Moved into position after routing of cables in overhead cable trays.

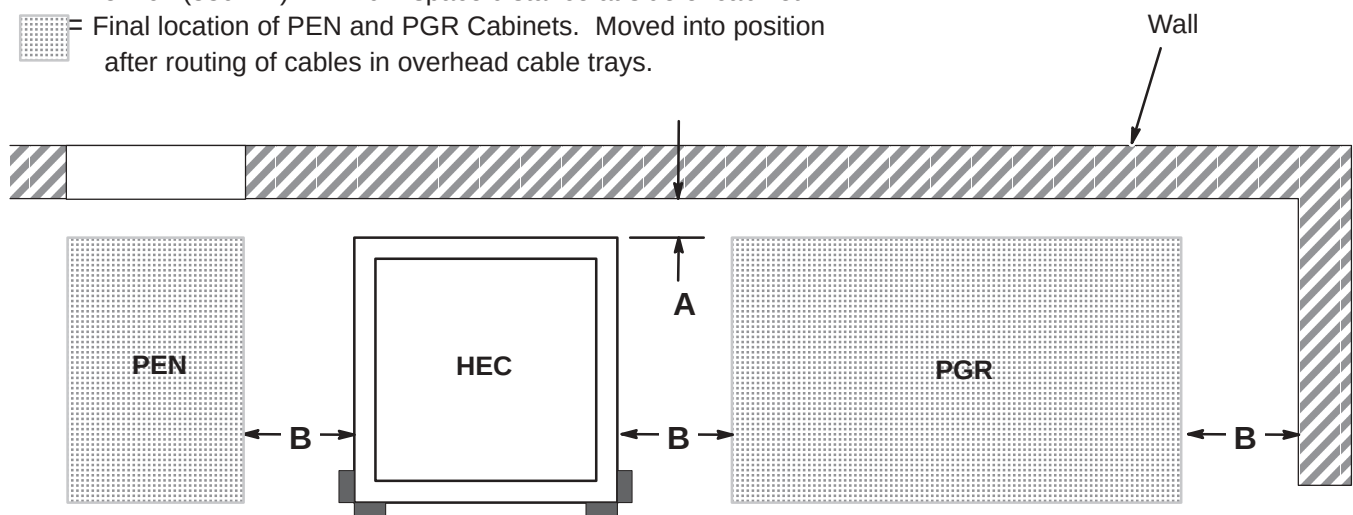2-2-2 CABINET WITH SHIPPING PALLET

The HEC cabinet is delivered to site secured to a pallet. The cabinet must be removed from the pallet and positioned in Equipment Room according to site layout plans. The cabinet moving dollies that were delivered to site on a pallet are required to move the HEC cabinet off the pallet.

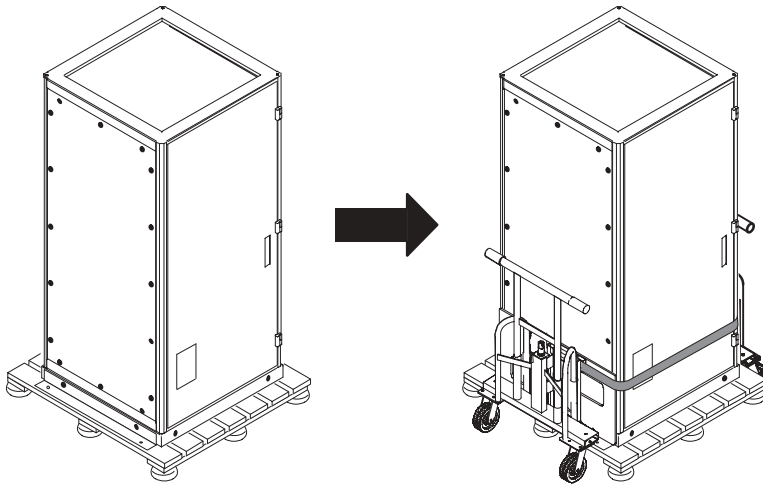
2-2-3 CABINET MOVING DOLLIES

The cabinet moving dollies have been shipped to site on a separate pallet. If dollies are not on site, they can be ordered online at www.umi-dollyshop.com. For technical questions, contact: UMI, 1-800-638-7457.

Note: The dollies have the following features to aid in their use:


- Wheel Locking Pin – Keeps the wheels locked into a certain direction
- Parking Brakes – For holding the cabinet in position
- Slide and Lock – The dollies can be removed in limited space (15 inches or 380mm)
- Twinning Hardware – hardware to connect each dolly to the other dolly during transport without a cabinet

1. Unfasten and remove dollies from shipping pallet.
2. If dollies need to be moved to different area where the HEC cabinet is located, attach the two dollies together using the twinning hardware. Refer to note above.

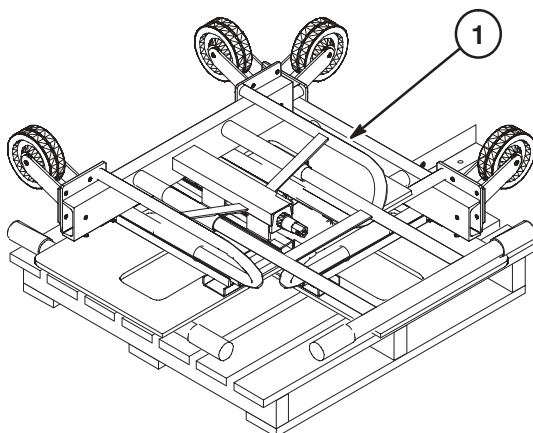
2-2-4 MOVE HEC CABINET INTO EQUIPMENT ROOM

1. Unfasten screws for each angle shipping bracket and remove bracket from each side of cabinet.
2. Install the six bushings, three on each side, to the bottom of the cabinet.

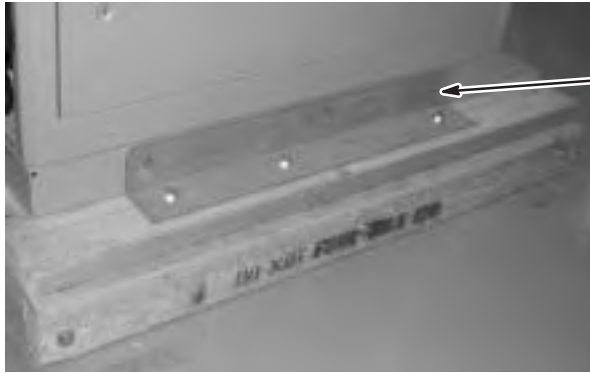
1

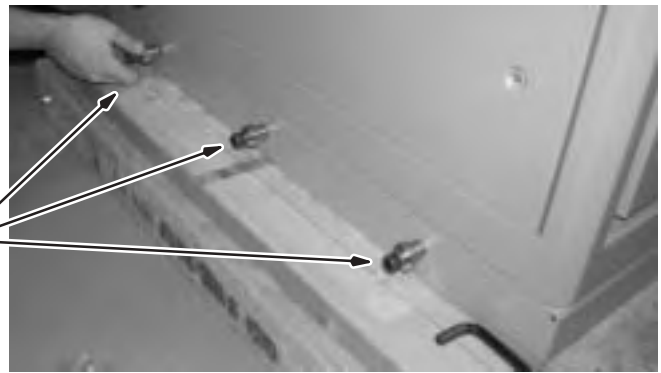
2

3. Maneuver the dolly so that all three bushings slip into the slide lock on the attachment plate.
4. Slide the dolly sideways to position the bushings in the lifting part of the attachment plate.

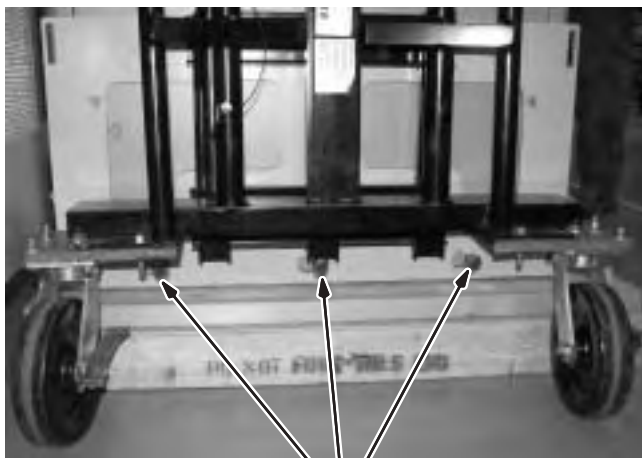
3

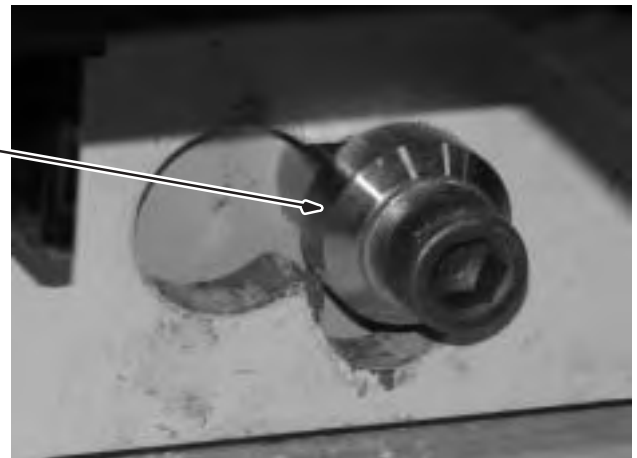
4

2-2-4 MOVE HEC CABINET INTO EQUIPMENT ROOM (continued)

5. Raise the dolly by turning the jack screw clockwise. This will lock the dolly into position.
6. Secure to dollies with strap. Make sure the strap goes through the slot of the attachment plate and the clamp is on flat surface of cabinet.
7. Tighten strap and close clamp.
8. Lift cabinet off of pallet by turning the jack screw clockwise until the cabinet surface is free from the pallet surface.
Note: Alternate raising each dolly in 1/4 inch (6mm) increments until the cabinet is above the pallet.

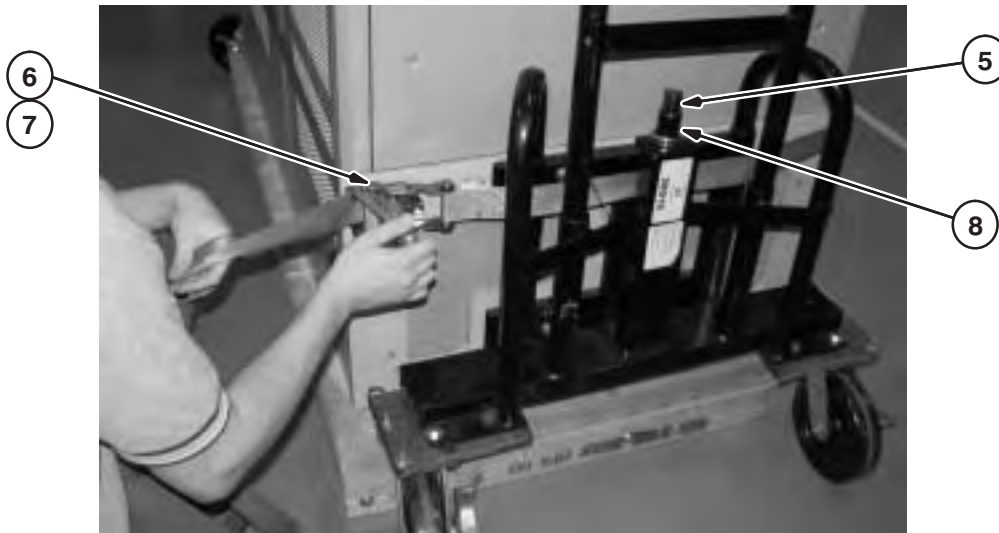
9. Slide the pallet out from under the cabinet.
10. Lower the cabinet by turning the jack screw counter-clockwise until the cabinet is about 1 inch (25mm) above the floor. **Note:** Alternate lowering each side of cabinet.
11. Move the cabinet into the equipment room.

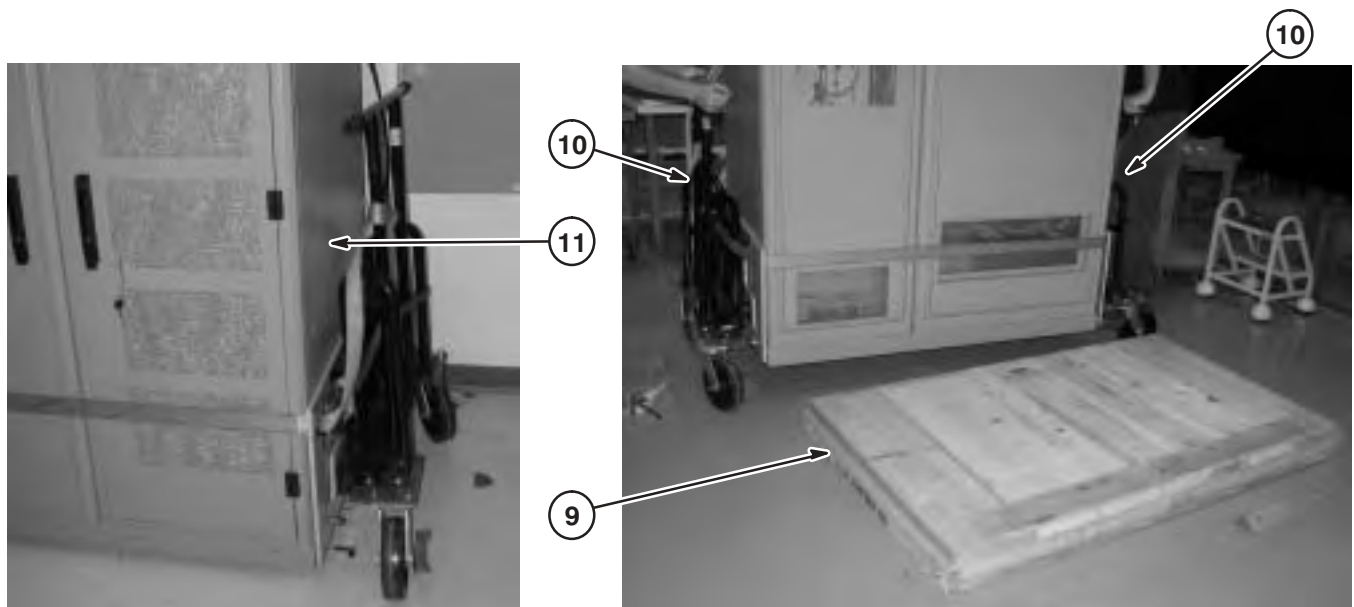2-2-4 MOVE HEC CABINET INTO EQUIPMENT ROOM (continued)

12. Refer to Section 2-2-1 for placement of cabinet and lower cabinet into position.
13. Leave dollies attached to cabinet. They will be used for moving and positioning the PGR and PEN cabinets when they are delivered to site.
14. Remove strap and store in bag delivered with dollies that contained the bushings.

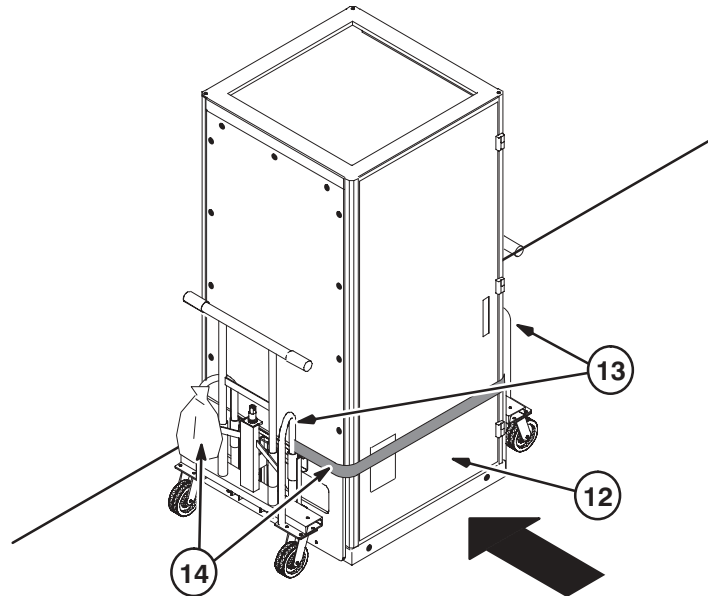2-3 INSTALLATION OF MAIN DISCONNECT PANEL

The MDP (Main Disconnect Panel, 5161287) is to be mounted and connected to facility power by customer electrician. Refer to architectural site layout drawings for specifications and mounting locations.

2-4 HEC TO CRYO-COOLER HOSE KIT

The HEC-CRY Hose Kit (5268715) is delivered to site and will be installed when Cryo-Cooler Compressor is delivered with magnet. Store the gray and yellow hoses (5264039 & 5264063), the Push-Lok fasteners (5272163 & 5268082), the Ball Valves (5272392), and the Flared Compression fittings (5268775) in safe place for installation later in the process.


2-5 HEC TO CRYO-COOLER HOSE INSULATION KIT

Along with the above hose kit, store the Flexible Foam Insulation Kit (5266250) in a safe place for until installed at a later time in the installation process.

2-6 CUSTOMER CONNECTIONS TO HEC CABINET

Customer is to provide and install liquid coolant, pipe/hoses, filters and connectors to the HEC cabinet.

Customer is also providing and installing facility power to the HEC cabinet from the MDP. Make sure to include a service loop for both the incoming power and the power cable routed to the Cryo-Cooler. The space above the power box has been provided for this service loop. Refer to site drawings for specifications.

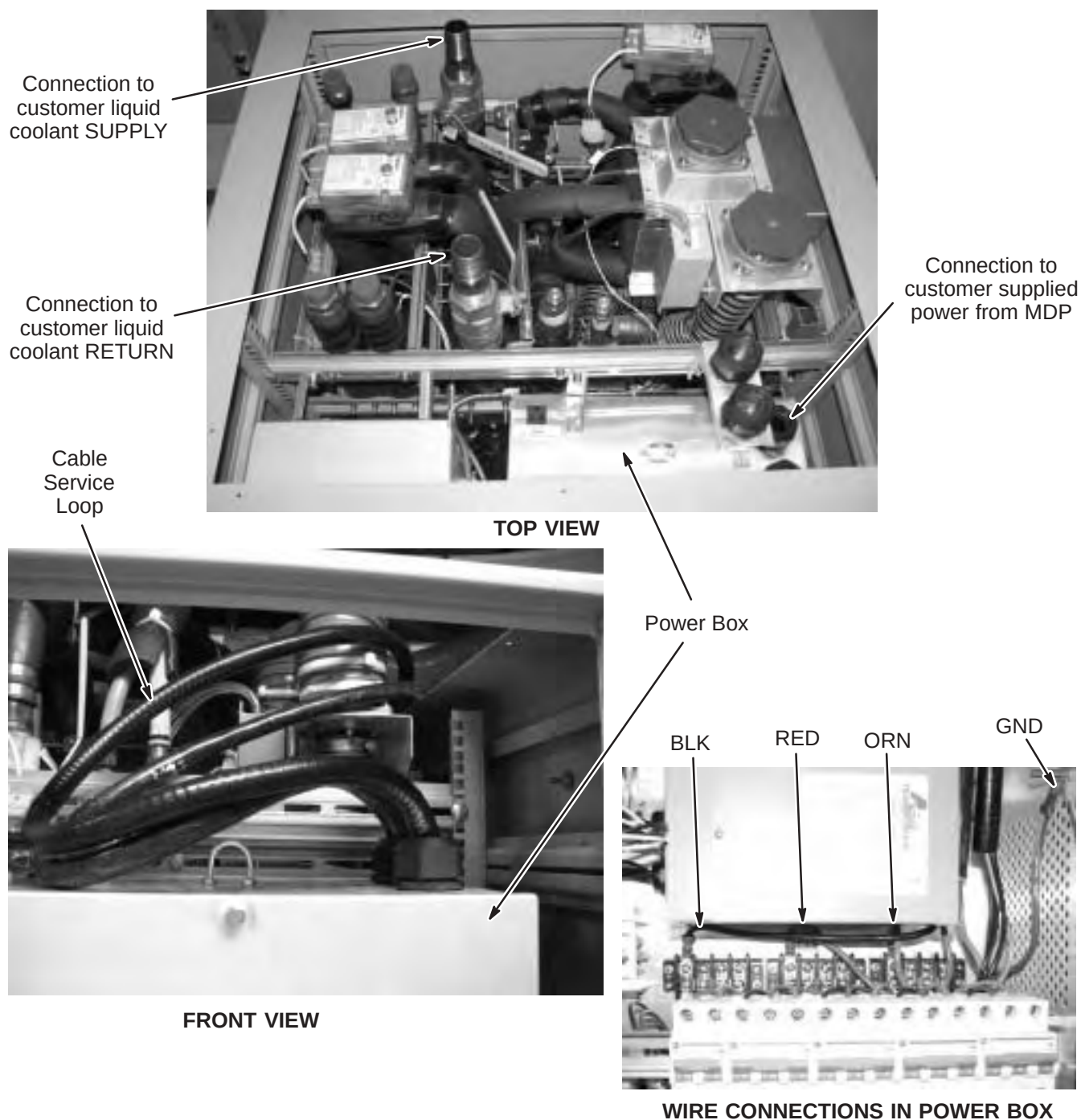2-7 INSTALL RF SCREEN ROOM WALL MOUNTING FRAMES AND SPW

Refer to architectural site layout drawings to determine location of SPW (Secondary Pen Wall) and Penetration Panel frames.


WARNING!

DUE TO WEIGHT OF SPW, AT LEAST TWO PEOPLE ARE REQUIRED TO LIFT AND MOUNT SPW INTO FRAME.

2-7-1 INSTALL MOUNTING FRAMES IN RF SCREEN ROOM OPENING

NOTE: The orientation of the Penetration Panel frame is different than the SPW frame.


1. From Equipment Room side of RF Screen Room wall, install the SPW and Penetration panel Frames. The Penetration Panel and SPW frames can be used as templates for hole locations.
IMPORTANT: Make sure frame is positioned so that gradient cable strain relief bracket mounting studs are at the top, left side.
2. Secure both frames to Scan Room wall. Customer to provide fasteners and instructions for installation to achieve RF integrity.

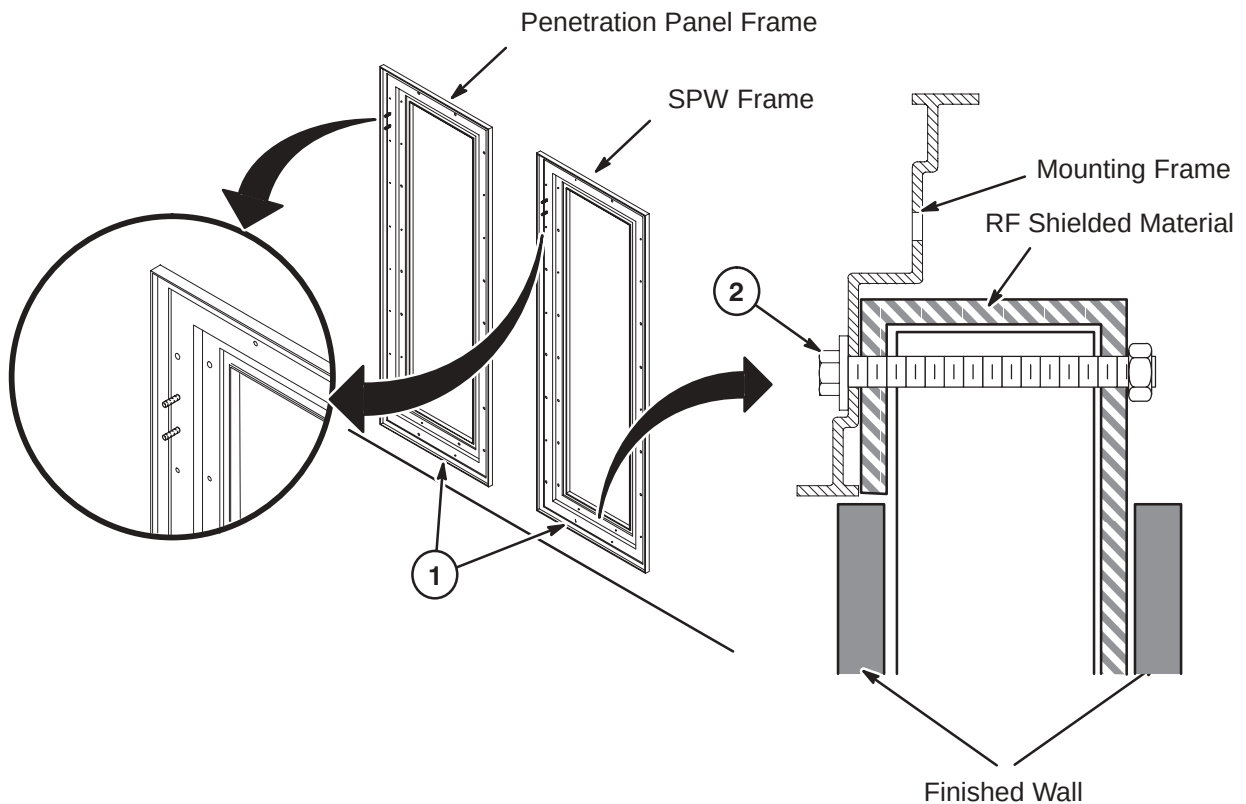2-7-2 INSTALL SPW

1. Move SPW into designated opening and mount to frame using 20 M5 x12mm SS Screws (2381274-2) and M5 Washers (2381270-6).
2. On Equipment Room side, install Gradient Cable Strain Relief Bracket (5182786) and secure on two studs using a M8 SS Hex Nut (2381089-6) on each stud.
3. Install three U-bolts (2113244) in bracket and loosely fasten with 3/8-16 Brass hex Nut (46-208935P10).
4. On Magnet Room side, place one 1 inch washer (46-251379P4) on each stud.
5. On Magnet Room side, install Gradient Cable Strain Relief Bracket (5182786) and secure on two studs using a M8 SS Hex Nut (2381089-6) on each stud.
6. Slide Air inlet Filter (5290917) on to end of J105 waveguide.
7. Secure Air Filter firmly in place with 4-inch Plastic Hose Clamp (5263772-3).


CAUTION: Do NOT slide filter mounting flange past end of waveguide. Damage to filter may result.

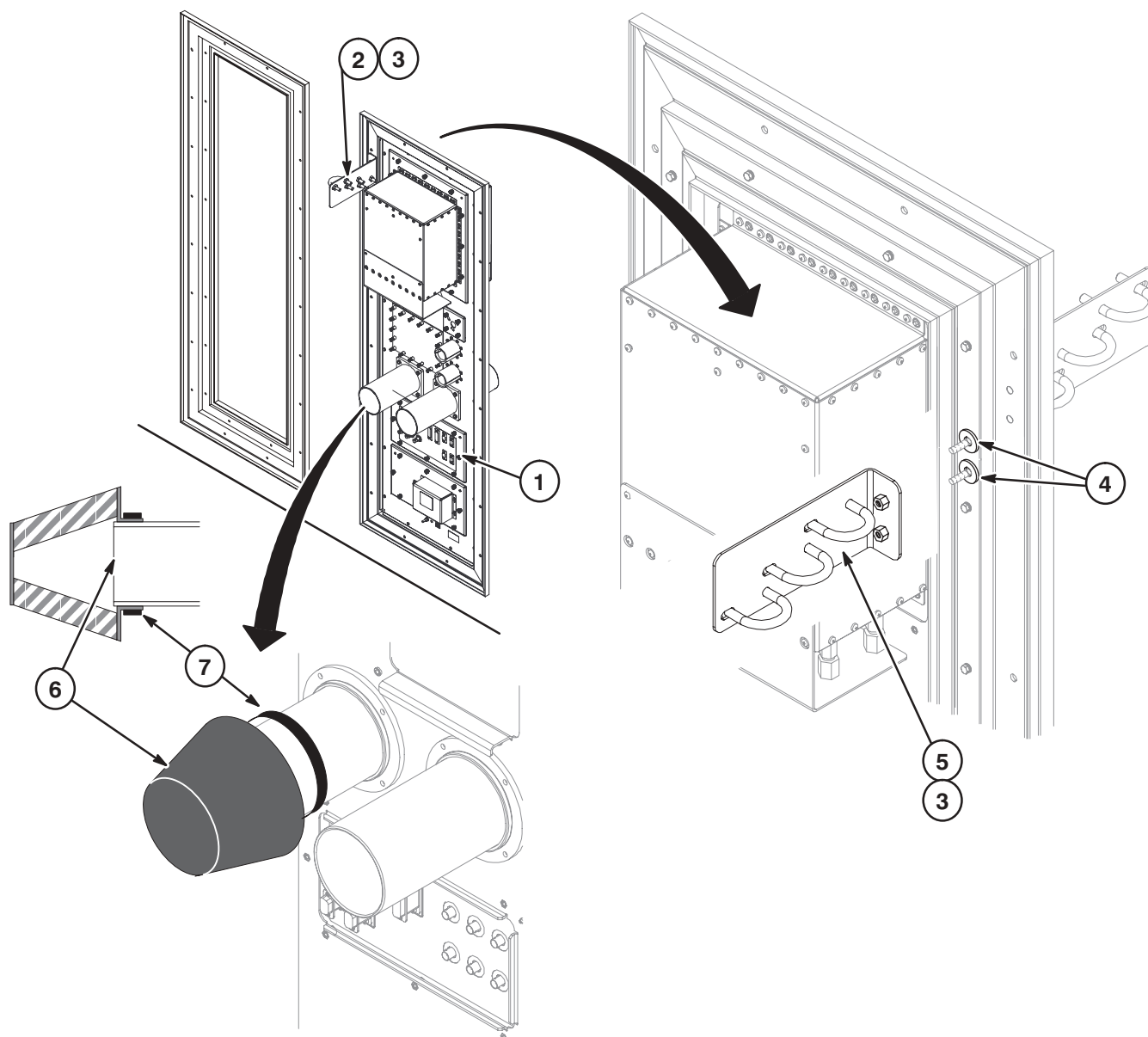2-8 INSTALL LOCK-OUT TAG-OUT PLACARD

The Lock-out Tag-out placard (5197122) must be secured to the wall next to the MDP.

1. Remove covering of double-sided tape on back side of placard.
2. Press and secure to surface next to MDP.

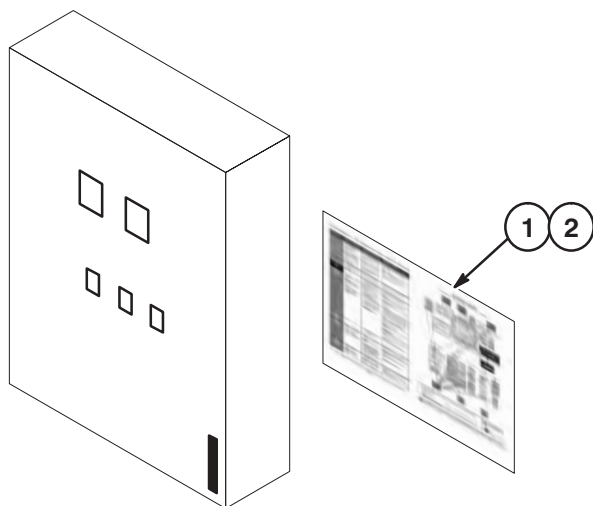
2-9 INSTALL MAGNET ENCLOSURE CABLE CURTAIN (Customer Option for MR750 only)

The Curtain Kit Collector (5241434) is delivered with this installation. It is provided if the customer site plans to install the curtain. This option is only installed for MR750.


If curtain is to be installed, the installation instructions are found in the following pages as part of 5267708DDW. The curtain kit is installed in two stages:

- First Stage – The Curtain Rail Front and Rear Assemblies (5237762 and 5237763) are to be installed prior to the delivery of the remaining system cabinets and cables, which will be delivered to site after magnet installation. Refer to sheet one of 5267708DDW on the next page for location dimensions of Curtain Rail from magnet iso-center or rear flange.
- Second Stage – The front and rear Curtain Covers (5241217 and 5241217-2) are installed after magnet enclosure has been installed and the Curtain Rail is in place in the ceiling. Refer to Direction 5500102, Section 5, for Curtain Cover installation procedures.


2-9 INSTALL MAGNET ENCLOSURE CABLE CURTAIN (continued)

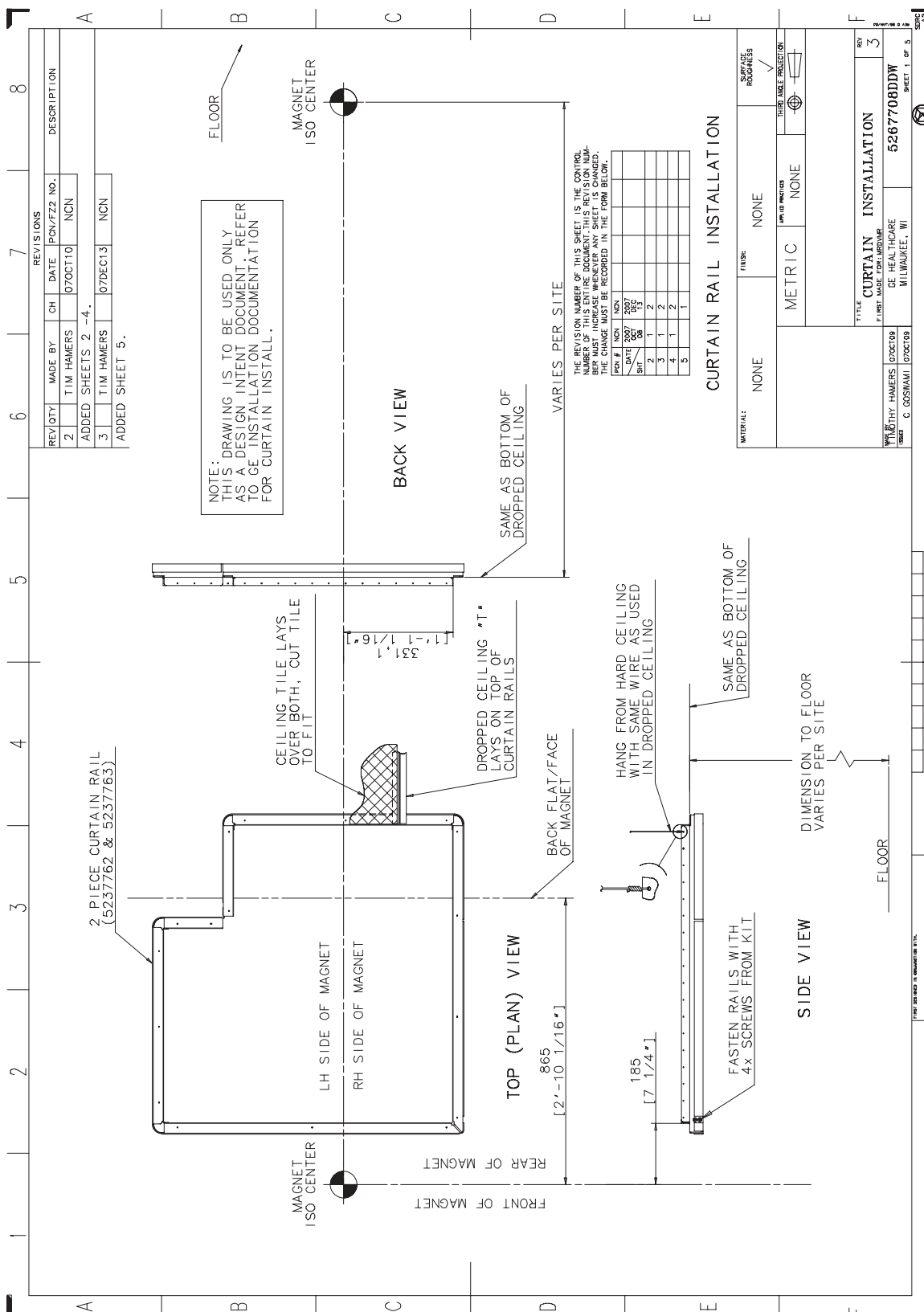
2-9 INSTALL MAGNET ENCLOSURE CABLE CURTAIN (continued)

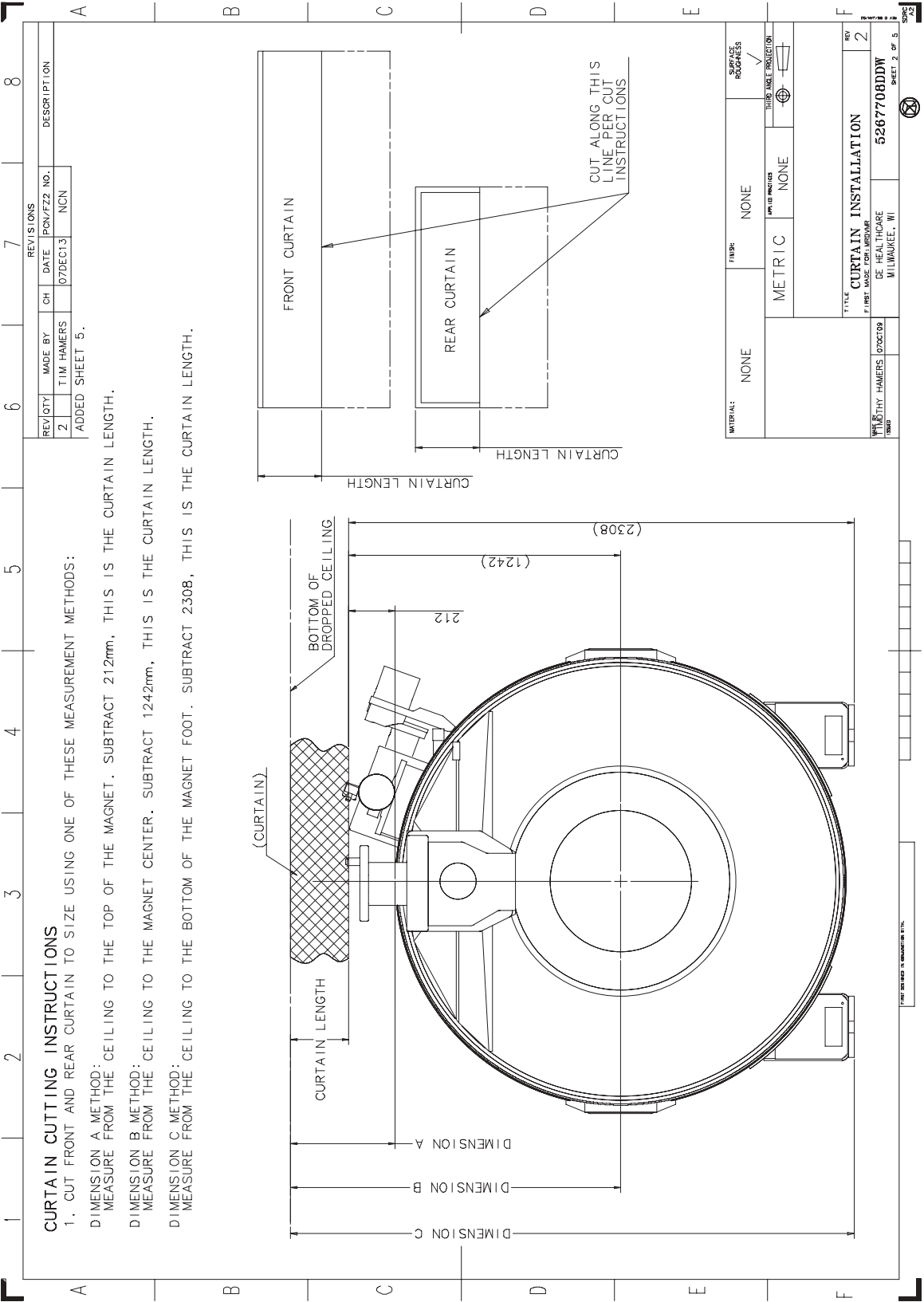
2-9 INSTALL MAGNET ENCLOSURE CABLE CURTAIN (continued)

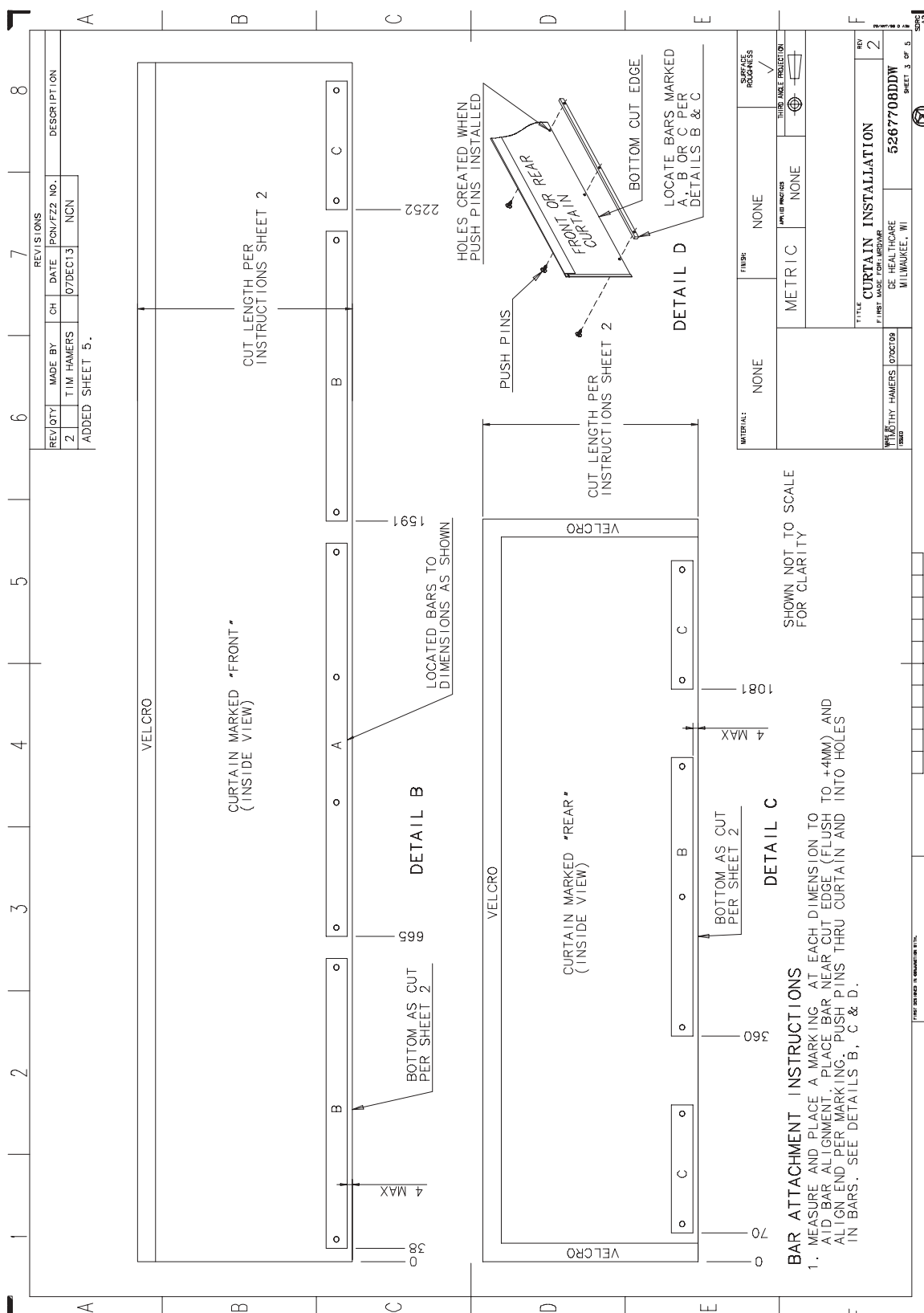2-9 INSTALL MAGNET ENCLOSURE CABLE CURTAIN (continued)

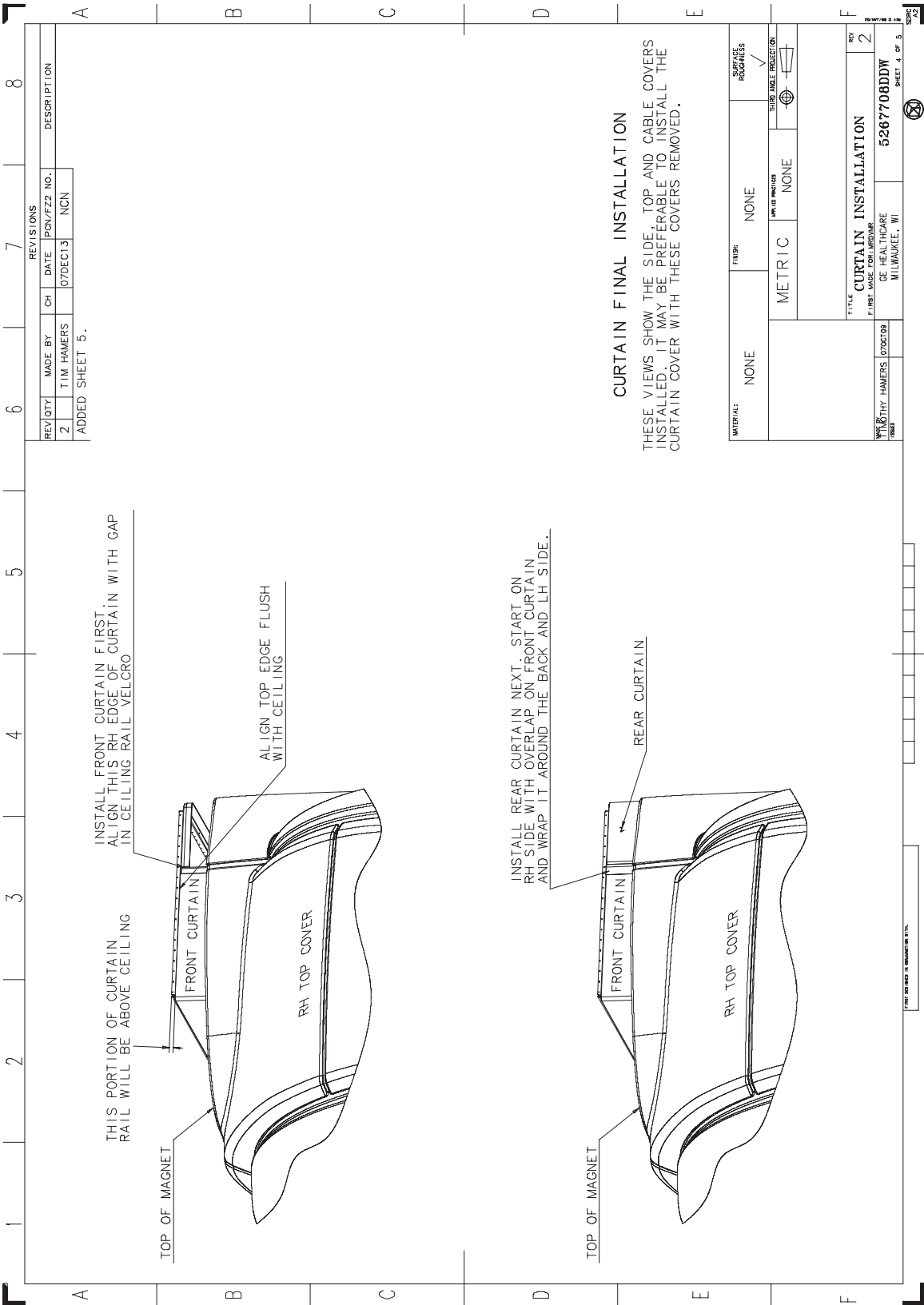2-9 INSTALL MAGNET ENCLOSURE CABLE CURTAIN (continued)

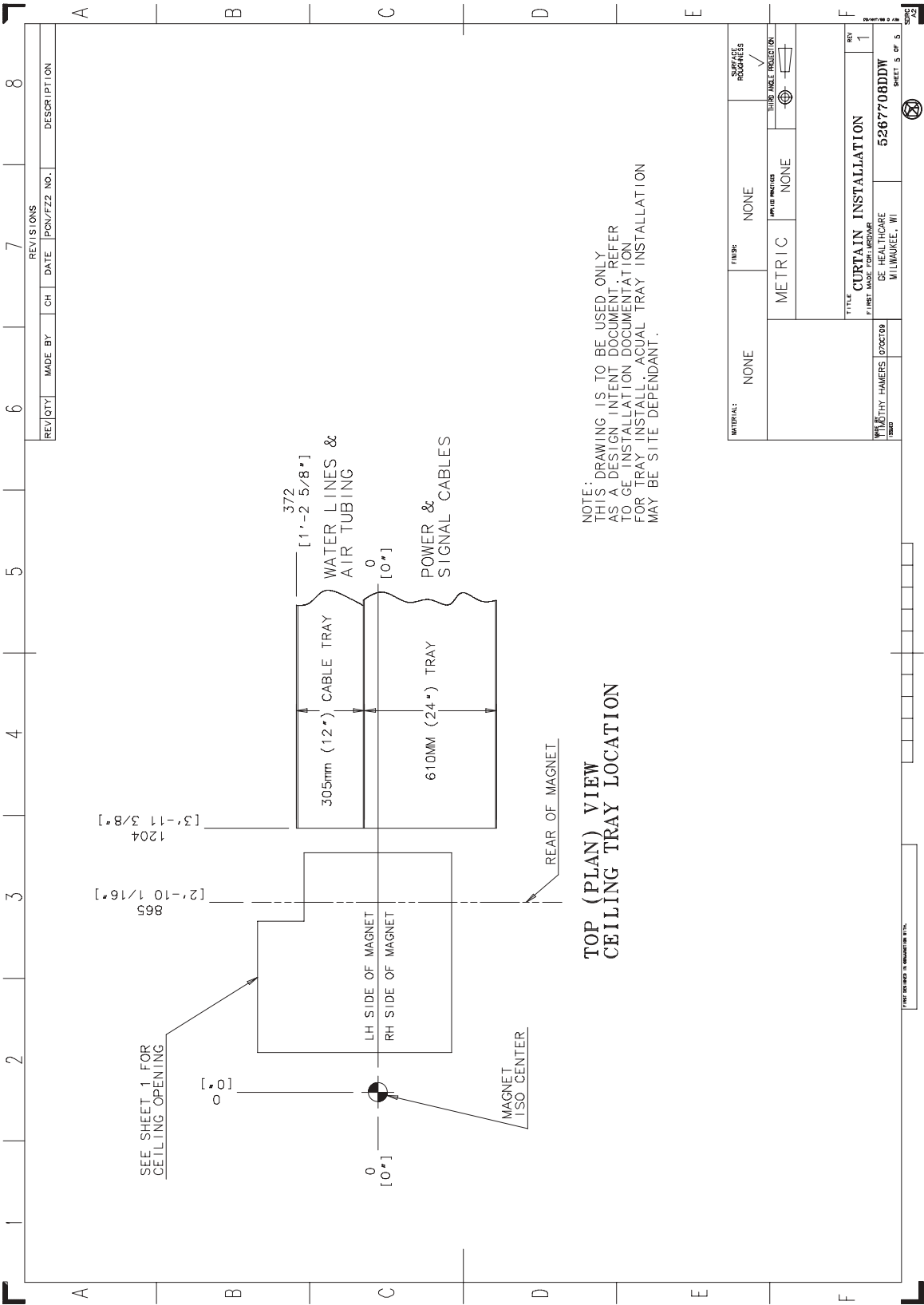SECTION 3 – COMPLETION CHECKLIST

TABLE OF CONTENTS

<u>SECTION</u>	<u>PAGE</u>
3-1 MECHANICAL INSTALLATION COMPLETION	3-2
3-2 GE FIELD ENGINEER RESPONSIBILITIES	3-2

3-1 MECHANICAL INSTALLATION COMPLETION

The following steps need to be reviewed and completed by GE Field Engineer, especially if the system was installed by non-GE personnel. Check off boxes at the front of each item below when complete.

- ☐ Check mechanical flowchart to make sure all steps are complete
- ☐ Resolve shipment shortages
- ☐ Resolve omissions by mechanical contractors
- ☐ Make sure all rating plates are installed
- ☐ Recheck all wiring connections using cable map
- ☐ Correct for wiring errors if necessary

3-2 GE FIELD ENGINEER RESPONSIBILITIES

1. Record and enter applicable data into applicable site configuration files and records.
2. Complete Product Locator information for all installed serialized components, new or updated, via either:
 - **(U.S. Only)** *FE Site Verification Web Site*
or
 - Process and return product locator installation cards for all serialized components to:
Product Locator File, P.O. Box 414, W-523, Milwaukee, WI 53201-0414

Refer to Section 1-8 for details on submitting Product Locator information.

Note

Failure to fill out and return Product Locator Cards may result in failure of your site to receive future FMI's.

3. Store the delivered site's set of service tools and spares kit in service cabinet at site.
4. Leave the site's set of Manuals and CD-ROM at site. Install Service and Operator documentation on the host computer. Organize and set up reference cabinet for the Manuals.
5. Locate any Material Safety Data Sheets (MSDS).. They must be retained on site. Customer is to be informed that material with MSDS was brought into site and customer should know/decide where MSDS should be retained at site.