



GE Medical Systems

Technical Publications

Direction 2397845
Revision 6

Signa® Software Option Installation With eLicense Generated Keys

Copyright® 2007 by General Electric Co.

Operating Documentation

Language Policy for Service Documentation

ПРЕДУПРЕЖДЕНИЕ (BG)	Това упътване за работа е налично само на английски език. - Ако доставчикът на услугата на клиента изиска друг език, задължение на клиента е да осигури превод. - Не използвайте оборудването, преди да сте се консултирали и разбрали упътването за работа. - Неспазването на това предупреждение може да доведе до нараняване на доставчика на услугата, оператора или пациента в резултат на токов удар, механична или друга опасност.
警告 (ZH-CN)	本维修手册仅提供英文版本。 - 如果客户的维修服务人员需要非英文版本，则客户需自行提供翻译服务。 - 未详细阅读和完全理解本维修手册之前，不得进行维修。 - 忽略本警告可能对维修服务人员、操作人员或患者造成电击、机械伤害或其他形式的伤害。
警告 (ZH-HK)	本服務手冊僅提供英文版本。 - 倘若客戶的服務供應商需要英文以外之服務手冊，客戶有責任提供翻譯服務。 - 除非已參閱本服務手冊及明白其內容，否則切勿嘗試維修設備。 - 不遵從本警告或會令服務供應商、網絡供應商或病人受到觸電、機械性或其他危險。
警告 (ZH-TW)	本維修手冊僅有英文版。 - 若客戶的維修廠商需要英文版以外的語言，應由客戶自行提供翻譯服務。 - 請勿試圖維修本設備，除非 您已查閱並瞭解本維修手冊。 - 若未留意本警告，可能導致維修廠商、操作員或病患因觸電、機械或其他危險而受傷。
UPOZORENJE (HR)	Ovaj servisni priručnik dostupan je na engleskom jeziku. - Ako davatelj usluge klijenta treba neki drugi jezik, klijent je dužan osigurati prijevod. - Ne pokušavajte servisirati opremu ako niste u potpunosti pročitali i razumjeli ovaj servisni priručnik. - Zanemarite li ovo upozorenje, može doći do ozljede davatelja usluge, operatera ili pacijenta uslijed strujnog udara, mehaničkih ili drugih rizika.
VÝSTRAHA (CS)	Tento provozní návod existuje pouze v anglickém jazyce. - V případě, že externí služba zákazníkům potřebuje návod v jiném jazyce, je zajištění překladu do odpovídajícího jazyka úkolem zákazníka. - Nesnažte se o údržbu tohoto zařízení, aniž byste si přečetli tento provozní návod a pochopili jeho obsah. - V případě nedodržování této výstrahy může dojít k poranění pracovníka prodejního servisu, obslužného personálu nebo pacientů vlivem elektrického proudu, respektive vlivem mechanických či jiných rizik.
ADVARSEL (DA)	Denne servicemanual findes kun på engelsk. - Hvis en kundes tekniker har brug for et andet sprog end engelsk, er det kundens ansvar at sørge for oversættelse. - Forsøg ikke at servicere udstyret uden at læse og forstå denne servicemanual. - Manglende overholdelse af denne advarsel kan medføre skade på grund af elektrisk stød, mekanisk eller anden fare for tekniker, operatøren eller patienten.
WAARSCHUWING (NL)	Deze onderhoudshandleiding is enkel in het Engels verkrijgbaar. - Als het onderhoudspersoneel een andere taal vereist, dan is de klant verantwoordelijk voor de vertaling ervan. - Probeer de apparatuur niet te onderhouden alvorens deze onderhoudshandleiding werd geraadpleegd en begrepen is. - Indien deze waarschuwing niet wordt opgevolgd, zou het onderhoudspersoneel, de operator of een patiënt gewond kunnen raken als gevolg van een elektrische schok, mechanische of andere gevaren.

WARNING (EN)	This service manual is available in English only. • If a customer's service provider requires a language other than english, it is the customer's responsibility to provide translation services. • Do not attempt to service the equipment unless this service manual has been consulted and is understood. • Failure to heed this warning may result in injury to the service provider, operator or patient from electric shock, mechanical or other hazards.
HOIATUS (ET)	See teenindusjuhend on saadaval ainult inglise keeles • Kui klienditeeninduse osutaja nõuab juhendit inglise keelest erinevas keeles, vastutab klient tõlketeenuse osutamise eest. • Ärge üritage seadmeid teenindada enne eelnevalt käesoleva teenindusjuhendiga tutvumist ja sellest aru saamist. • Käesoleva hoiatuse eiramine võib põhjustada teenuseosutaja, operaatori või patsiendi vigastamist elektrilöögi, mehaanilise või muu ohu tagajärjel.
VAROITUS (FI)	Tämä huolto-ohje on saatavilla vain englanniksi. • Jos asiakkaan huoltohenkilöstö vaatii muuta kuin englanninkielistä materiaalia, tarvittavan käännöksen hankkiminen on asiakkaan vastuulla. • Älä yritä korjata laitteistoa ennen kuin olet varmasti lukenut ja ymmärtänyt tämän huolto-ohjeen. • Mikäli tätä varoitusta ei noudateta, seurauksena voi olla huoltohenkilöstön, laitteiston käyttäjän tai potilaan vahingoittuminen sähköiskun, mekaanisen vian tai muun vaaratilanteen vuoksi.
ATTENTION (FR)	Ce manuel d'installation et de maintenance est disponible uniquement en anglais. • Si le technicien d'un client a besoin de ce manuel dans une langue autre que l'anglais, il incombe au client de le faire traduire. • Ne pas tenter d'intervenir sur les équipements tant que ce manuel d'installation et de maintenance n'a pas été consulté et compris. • Le non-respect de cet avertissement peut entraîner chez le technicien, l'opérateur ou le patient des blessures dues à des dangers électriques, mécaniques ou autres.
WARNUNG (DE)	Diese Serviceanleitung existiert nur in englischer Sprache. • Falls ein fremder Kundendienst eine andere Sprache benötigt, ist es Aufgabe des Kunden für eine entsprechende Übersetzung zu sorgen. • Versuchen Sie nicht diese Anlage zu warten, ohne diese Serviceanleitung gelesen und verstanden zu haben. • Wird diese Warnung nicht beachtet, so kann es zu Verletzungen des Kundendiensttechnikers, des Bedieners oder des Patienten durch Stromschläge, mechanische oder sonstige Gefahren kommen.
ΠΡΟΕΙΔΟΠΟΙΗΣΗ (EL)	Το παρόν εγχειρίδιο σέρβις διατίθεται μόνο στα αγγλικά. • Εάν ο τεχνικός σέρβις ενός πελάτη απαιτεί το παρόν εγχειρίδιο σε γλώσσα εκτός των αγγλικών, αποτελεί ευθύνη του πελάτη να παρέχει τις υπηρεσίες μετάφρασης. • Μην επιχειρήσετε την εκτέλεση εργασιών σέρβις στον εξοπλισμό αν δεν έχετε συμβουλευτεί και κατανοήσει το παρόν εγχειρίδιο σέρβις. • Αν δεν προσέξετε την προειδοποίηση αυτή, ενδέχεται να προκληθεί τραυματισμός στον τεχνικό σέρβις, στο χειριστή ή στον ασθενή από ηλεκτροπληξία, μηχανικούς ή άλλους κινδύνους.
FIGYELMEZTETÉS (HU)	Ezen karbantartási kézikönyv kizárólag angol nyelven érhető el. • Ha a vevő szolgáltatója angoltól eltérő nyelvre tart igényt, akkor a vevő felelőssége a fordítás elkészítése. • Ne próbálja elkezdni használni a berendezést, amíg a karbantartási kézikönyvben leírtakat nem értelmezték. • Ezen figyelmeztetés figyelmen kívül hagyása a szolgáltató, működtető vagy a beteg áramütés, mechanikai vagy egyéb veszélyhelyzet miatti sérülését eredményezheti.
AÐVÖRUN (IS)	Þessi þjónustuhandbók er aðeins ááenleg á ensku. • Ef að þjónustuveitandi viðskiptamanns þarfnast annas tungumáls en ensku, er það skylda viðskiptamanns að skaffa tungumálþjónustu. • Reynið ekki að afgreiða tækið nema að þessi þjónustuhandbók hefur verið skoðuð og skilin. • Brot á sinna þessari aðvörun getur leitt til meiðsla á þjónustuveitanda, stjórnaða eða sjúklings frá raflosti, vélrænu eða öðrum áhættum.

AVVERTENZA
(IT)

Il presente manuale di manutenzione è disponibile soltanto in lingua inglese.

- Se un addetto alla manutenzione richiede il manuale in una lingua diversa, il cliente è tenuto a provvedere direttamente alla traduzione.
- Procedere alla manutenzione dell'apparecchiatura solo dopo aver consultato il presente manuale ed averne compreso il contenuto.
- Il mancato rispetto della presente avvertenza potrebbe causare lesioni all'addetto alla manutenzione, all'operatore o ai pazienti provocate da scosse elettriche, urti meccanici o altri rischi.

警告
(JA)

このサービスマニュアルには英語版しかありません。

- サービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。
- このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないでください。
- この警告に従わない場合、サービスを担当される方、操作員あるいは患者

さんが、感電や機械的又はその他の危険により負傷する可能性があります。

경고
(KO)

본 서비스 매뉴얼은 영어로만 이용하실 수 있습니다.

- 고객님의 서비스 제공자가 영어 이외의 언어를 요구할 경우, 번역 서비스를 제공하는 것은 고객님의 책임입니다.
- 본 서비스 매뉴얼을 참조하여 숙지하지 않은 이상 해당 장비를 수리하려고 시도하지 마십시오.
- 본 경고 사항에 유의하지 않으면 전기 쇼크, 기계적 위험, 또는 기타 위험으로 인해 서비스 제공자, 사용자 또는 환자에게 부상을 입힐 수 있습니다.

BRĪDINĀJUMS
(LV)

Šī apkopes rokasgrāmata ir pieejama tikai anglū valodā.

- Ja klienta apkopes sniedzējam nepieciešama informācija citā valodā, klienta pienākums ir nodrošināt tulkojumu.
- Neveiciet aprīkojuma apkopi bez apkopes rokasgrāmatas izlasīšanas un saprašanas.
- Šī brīdinājuma neievērošanas rezultātā var rasties elektriskās strāvas trieciena, mehānisku vai citu faktoru izraisītu traumu risks apkopes sniedzējam, operatoram vai pacientam.

[SPĖJIMAS
(LT)

Šis eksploatavimo vadovas yra tik anglų kalba.

- Jei kliento paslaugų tiekėjas reikalauja vadovo kita kalba – ne anglų, suteikti vertimo paslaugas privalo klientas.
- Nemėginkite atlikti įrangos techninės priežiūros, jei neperskaitėte ar nesupratote šio eksploatavimo vadovo.
- Jei nepaisysite šio įspėjimo, galimi paslaugų tiekėjo, operatoriaus ar paciento sužalojimai dėl elektros šoko, mechaninių ar kitų pavojų.

ADVARSEL
(NO)

Denne servicehåndboken finnes bare på engelsk.

- Hvis kundens serviceleverandør har bruk for et annet språk, er det kundens ansvar å sørge for oversettelse.
- Ikke forsøk å reparere utstyret uten at denne servicehåndboken er lest og forstått.
- Manglende hensyn til denne advarselen kan føre til at serviceleverandøren, operatøren eller pasienten skades på grunn av elektrisk støt, mekaniske eller andre farer.

OSTRZEŻENIE
(PL)

Niniejszy podręcznik serwisowy dostępny jest jedynie w języku angielskim.

- Jeśli serwisant klienta wymaga języka innego niż angielski, zapewnienie usługi tłumaczenia jest obowiązkiem klienta.
- Nie próbować serwisować urządzenia bez zapoznania się z niniejszym podręcznikiem serwisowym i zrozumienia go.
- Niezastosowanie się do tego ostrzeżenia może doprowadzić do obrażeń serwisanta, operatora lub pacjenta w wyniku porażenia prądem elektrycznym, zagrożenia mechanicznego bądź innego.

AVISO (PT-BR)	Este manual de assistência técnica encontra-se disponível unicamente em inglês. • Se outro serviço de assistência técnica solicitar a tradução deste manual, caberá ao cliente fornecer os serviços de tradução. • Não tente reparar o equipamento sem ter consultado e compreendido este manual de assistência técnica. • A não observância deste aviso pode ocasionar ferimentos no técnico, operador ou paciente decorrentes de choques elétricos, mecânicos ou outros.
ATENÇÃO (PT-PT)	Este manual de assistência técnica só se encontra disponível em inglês. • Se qualquer outro serviço de assistência técnica solicitar este manual noutro idioma, é da responsabilidade do cliente fornecer os serviços de tradução. • Não tente reparar o equipamento sem ter consultado e compreendido este manual de assistência técnica. • O não cumprimento deste aviso pode colocar em perigo a segurança do técnico, do operador ou do paciente devido a choques elétricos, mecânicos ou outros.
ATENȚIE (RO)	Acest manual de service este disponibil doar în limba engleză. • Dacă un furnizor de servicii pentru clienți necesită o altă limbă decât cea engleză, este de datoria clientului să furnizeze o traducere. • Nu încercați să reparați echipamentul decât ulterior consultării și înțelegerii acestui manual de service. • Ignorarea acestui avertisment ar putea duce la rănirea depanatorului, operatorului sau pacientului în urma pericolelor de electrocutare, mecanice sau de altă natură.
ОСТОРОЖНО! (RU)	Данное руководство по техническому обслуживанию представлено только на английском языке. • Если сервисному персоналу клиента необходимо руководство не на английском, а на каком-то другом языке, клиенту следует самостоятельно обеспечить перевод. • Перед техническим обслуживанием оборудования обязательно обратитесь к данному руководству и поймите изложенные в нем сведения. • Несоблюдение требований данного предупреждения может привести к тому, что специалист по техобслуживанию, оператор или пациент получит удар электрическим током, механическую травму или другое повреждение.
UPOZORENJE (SR)	Ovo servisno uputstvo je dostupno samo na engleskom jeziku. • Ako klijentov serviser zahteva neki drugi jezik, klijent je dužan da obezbedi prevodilačke usluge. • Ne pokušavajte da opravite uređaj ako niste pročitali i razumeli ovo servisno uputstvo. • Zanemarivanje ovog upozorenja može dovesti do povređivanja serviser, rukovaoca ili pacijenta usled strujnog udara ili mehaničkih i drugih opasnosti.
UPOZORNENIE (SK)	Tento návod na obsluhu je k dispozícii len v angličtine. • Ak zákazníkovi poskytovateľ služieb vyžaduje iný jazyk ako angličtinu, poskytnutie prekladateľských služieb je zodpovednosťou zákazníka. • Nepokúšajte sa o obsluhu zariadenia, kým si neprečítate návod na obsluhu a neporozumiete mu. • Zanedbanie tohto upozornenia môže spôsobiť zranenie poskytovateľa služieb, obsluhujúcej osoby alebo pacienta elektrickým prúdom, mechanické alebo iné ohrozenie.
ATENCION (ES)	Este manual de servicio sólo existe en inglés. • Si el encargado de mantenimiento de un cliente necesita un idioma que no sea el inglés, el cliente deberá encargarse de la traducción del manual. • No se deberá dar servicio técnico al equipo, sin haber consultado y comprendido este manual de servicio. • La no observancia del presente aviso puede dar lugar a que el proveedor de servicios, el operador o el paciente sufran lesiones provocadas por causas eléctricas, mecánicas o de otra naturaleza.

VARNING
(SV)

Den här servicehandboken finns bara tillgänglig på engelska. .

- Om en kunds servicetekniker har behov av ett annat språk än engelska, ansvarar kunden för att tillhandahålla översättningstjänster.
- Försök inte utföra service på utrustningen om du inte har läst och förstår den här servicehandboken.
- Om du inte tar hänsyn till den här varningen kan det resultera i skador på serviceteknikern, operatören eller patienten till följd av elektriska stötar, mekaniska faror eller andra faror.

DİKKAT
(TR)

Bu servis kılavuzunun sadece ingilizcesi mevcuttur.

- Eğer müşteri teknisyeni bu kılavuzu ingilizce dışında bir başka lisandan talep ederse, bunu tercüme ettirmek müşteriye düşer.
- Servis kılavuzunu okuyup anlamadan ekipmanlara müdahale etmeyiniz.
- Bu uyarıya uyulmaması, elektrik, mekanik veya diğer tehlikelerden dolayı teknisyen, operatör veya hastanın yaralanmasına yol açabilir.

REVISION HISTORY

REV	DATE	PRIMARY REASON FOR CHANGE
0.....	Oct 28 2003	Initial version for Release.
1.....	Jun 30 2004	Addition of 3T, HFO4, and MFO4 options, update to Excite II options, addition of copy to CD method and site selection screen
2.....	Feb 10 2006	Added HD and HDx options, added model type Table 2-1
3.....	March 3, 2006	Added Important statement in section 3 to make sure install options matches the upgrade being performed.
4.....	July 19, 2006	Added section detailing how to reconfigure eLicense on systems that have been upgraded to HDx. Updated the Downloading Licenses section (now 2-3-1). Updated 3-2 to reflect the burning options available with Windows XP computers. Minor grammatical and style changes throughout.
5...	November 20, 2007	Added DVMR options and model type.
6.....	October 10, 2009	Added MR750, MR450, HDxt, and HDi options and model type to Table 2-1. Update Language Policy.

PAGE	REV	PAGE	REV	PAGE	REV	PAGE	REV
Title Page.....	2	Title Page	3	Title Page	4	Title Page.....	5
2-29.....	2	3-17	3	2-4	4	10.....	5
		3-20	3			23-42.....	5

PAGE	REV
Title Page.....	6
2-6.....	6
9 and 17	6

* This revision/letter corresponds to the indicated document's revision control system.

TABLE OF CONTENTS


TABLE OF CONTENTS	8
1- Overview.....	9
2- elicense Website Tasks	10
2-1 Create Licenses.....	10
2-1-1 Obtaining the System's Host ID	10
2-1-2 Adding a System to an Order (For New Installations)	11
2-1-3 Adding a System to an Order (For Upgrades).....	15
2-2 Viewing existing licenses for a system	17
2-3 Distributing Licenses.....	18
2-3-1 Downloading Licenses.....	19
2-3-2 Emailing Licenses	20
2-3-3 Saving Licenses to a Disk.....	20
2-3-4 Printing Licenses.....	22
2-4 Modifying System Information	22
3- Option Installation.....	24
3-1 eLicense Floppy.....	25
3-2 eLicense CD-ROM.....	28
3-2-1 Making the eLicense CD-ROM	28
3-2-2 Loading options from eCD	29
3-3 Option Key Auto Download	30
3-4 FTP Option keys to scanner	31
3-5 Manual Option Key Entry.....	31
4- Option Key Catalog Matrix	33
5- Glossary	50

1- OVERVIEW

This procedure describes how to generate and install software options on the Signa scanners running 11.x software or greater on Core products and OpenMR4 or greater software on Open products using eLicense. eLicense allows the field to generate license keys for MR scanners and removes the need to send out physical media for all options on the order. The first step in creating option keys is to go to the eLicense website and register the scanner. In order to access the eLicense website, you must have a computer with connectivity, as well as a single sign-on account. You can access the eLicense website at <http://elicense.gehealthcare.com>.

The MR system must be tied to the global order number (GON) through the host ID and system ID of the scanner. Once the site information is entered, the keys can be generated and distributed by: email, downloading, saving to an eFloppy or eCD, or printing. This file can then be installed on the system, through inSite download, eFloppy, eCD, or manual entry.

Illustration 1-1 shows the eLicense process flow. For additional information on how to generate and obtain these keys from the eLicense server, please refer to the “**eLicense User Guide**” under the **HELP** tab of the eLicense website.

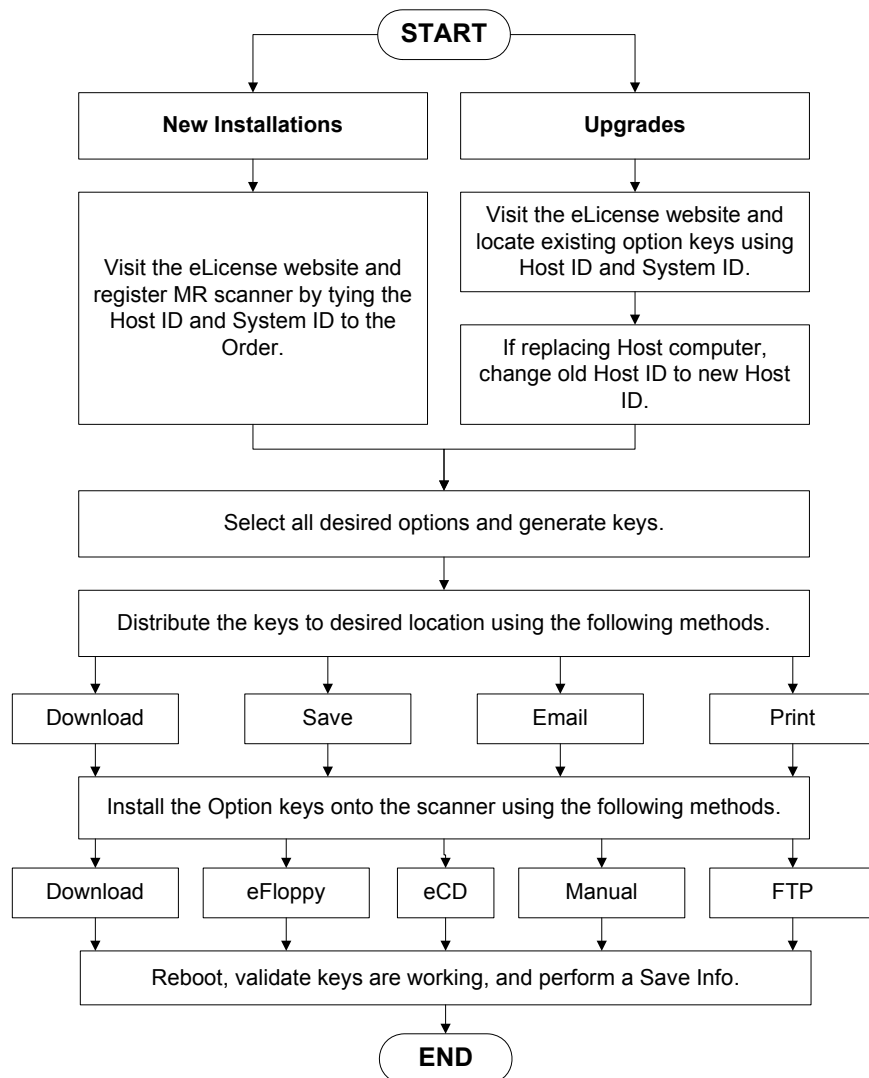

PROCESS FLOW OF OPTION KEY GENERATION
ILLUSTRATION 1-1


2- eLICENSE WEBSITE TASKS

2-1 Create Licenses

This feature creates option keys for a new or existing system based on an order. To create new licenses, you must have the order number(s) that lists the options ordered, as well as the host ID and system ID of the system.

2-1-1 Obtaining the System's Host ID

1. Select the **[Service Tools Desktop]**.
2. Select **[Guided Install]** from the Service Desktop Manager and press **[Start...]** to open Guided Install.
3. A command window displays the message **"Type Root Password"**. Type the system's root password and then press **<Enter>**. The Signa default root password is **"operator"**. If the customer has changed the root password, you need to obtain the correct root password from the customer.
4. Once the Guided Install window opens, select the *Option Keys* tab. The host ID is in the *Install Temporary Service Key* section. See Illustration 2-1.

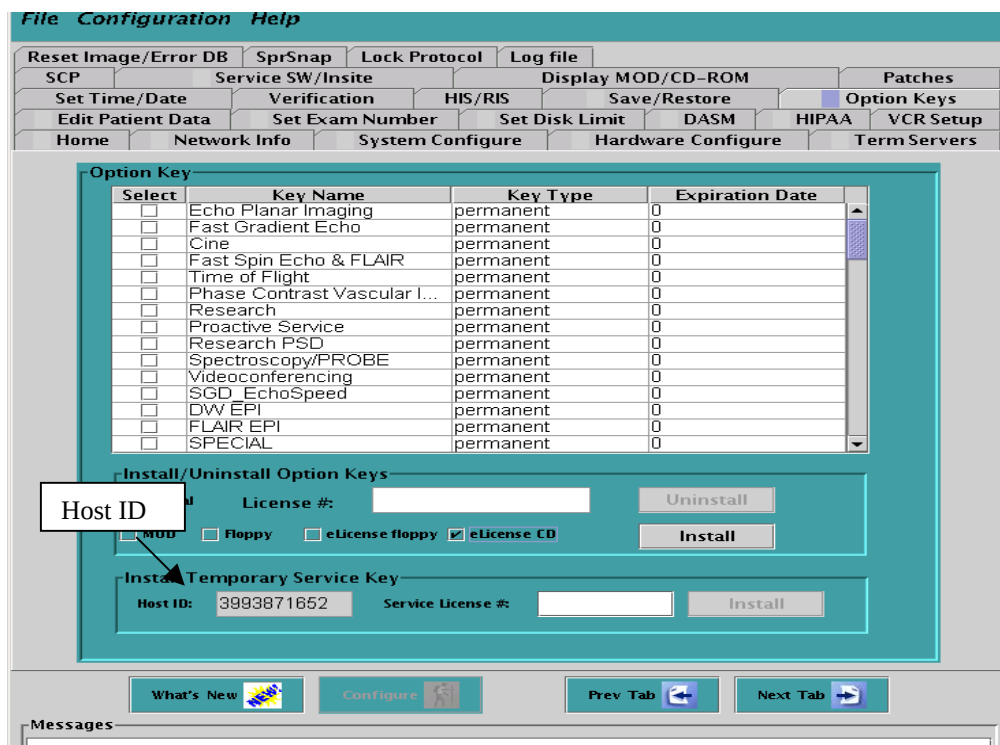

GUIDED INSTALL WINDOW WITH HOST ID
ILLUSTRATION 2-1

5. Another method of obtaining the host ID is to type **"sysinfo"** at a command prompt.

Note

"hostid" is not the valid command to get your Host ID. This causes an error when trying to generate keys, due to the incorrect host ID. You **MUST** use the methods listed above to get a valid ID.

2-1-2 Adding a System to an Order (For New Installations)

1. Go to the eLicense website <http://elicense.gehealthcare.com> and log on using your single sign-on username and password.
2. From the eLicense home page, enter the order number into the **Order Number** field and click **<Go>**. See Illustration 2-2. The order associated with the number entered from the eOM system displays. The order number may be a new global order number, or an FDO number, depending on when the order was placed. This number is listed on the packing list included with the order.

Address <http://apps.gemedicalsystems.com/elicense/index.jsp>

GE Medical Systems eLicense

Home | Admin | User Guide | Feedback Tuesday | October 7, 2003

Welcome Sadie

Order Management

Please enter an Order Number. If your order number includes a space or a dash, please include that as well.

Order Number

Examples:
1234-567890 or 1234 567890

Optional:
Enter a Host ID, System ID, or Serial Number and choose a modality to automatically add a system to the order.

Host ID

Choose Modality

System Configuration

Please enter a host id, system id, or serial number and choose a modality to view a list of a system's generated licenses.

Host ID

Choose Modality

*Note format options:
FDO Number or
new Oracle order #*

ELICENSE HOME PAGE
ILLUSTRATION 2-2

3. It is possible that the order number is not unique to your customer. If this is the case, a list of all matching sites displays. Illustration 2-3 shows this Site Selection screen. Select the desired system and click **<Go>**. This displays the order (See Illustration 2-4). If the order number is unique, eLicense automatically brings up the order and displays the *Order Information* screen.

Order Information

Go

☐ **Customer Name:** GE YOKOGAWA MEDICAL SYSTEM LTD
Ship-to Address: DAB在庫, 日野市旭が丘 4 - 7 - 1 2 7
City, State, Zip: HINO-SHI,,191-8503
Country: JP
Order Source: Internal

☒ **Customer Name:** WINTHROP UNIV HOSPITAL
Ship-to Address: SECOND ST, RECEIVING
City, State, Zip: MINEOLA, NY, 11501
Country: US
Order Source: EOM

Go

SITE SELECTION SCREEN
ILLUSTRATION 2-3

- a. Verify that the Customer Information at the top left of the screen is the correct customer.
4. Click **<Add a System>** to associate the order to a specific system. See Illustration 2-4.

Order Information

FDO: 0002 002612 /
Site: CHILDRENS HOSPITALS & CLINICS
Billing Number: 000146914
Address: ACCOUNTS PAYABLE
Street: 2910 CENTRE POINT DR
City, State Zip: ROSEVILLE MN 55113

Add a System

1. Verify Customer Information

2. Add a System Tab

*Note Order Source:
Your order of interest will
most likely be an "EOM"
or "Online" Source*

System Number	Description	Quantity	Unit Price	Total Price
M3033KK	APA2 FOR FORWARD PR	1	0	0
M3090KJ	2D FIESTA FOR EXCH	1	0	0
S7502JS	CONNECTPROPLUSBARC	1	0	0
M3090PL	CONNECTPRO FOR EX	1	0	0
M1090PM	BAR CODE READER/CON	1	0	0
M1033PA	PERFORMED PROCEDURE	1	0	0
M3033MV	IDRIVE PRO PLUS FOR EXCH	1	0	0
M80501DA	OM-AW4.1-LINUX-USA	1	0	0
S7501AF	DUAL ARRAY	1	0	0
M3087N	NON OBSOLESCENCE	1	0	0
M1000NT	18" COLOR LCD MONITOR	1	0	0
M1087SP	CTL COIL	1	0	0
S7502NB	8CH NV NON-OB FULFILL	1	0	0
M1000MN	KEYBOARD	1	0	0
M1000MV	TABLE	1	0	0
M1085KC	CHILLER FOR WIDE OPEN 60H	1	0	0
M1000NV	GLOBAL MODEM KIT	1	0	0
M1080M	SUMTOMO COMPRESSOR KIT	1	0	0

ORDER INFORMATION PAGE
ILLUSTRATION 2-4

5. Enter the Host ID obtained in step 2-1-1, select Host ID from the menu, select the MR modality, and then click **<Go>**. See Illustration 2-5.

Enter System Information

Host ID  

ENTER SYSTEM INFORMATION POP-UP
ILLUSTRATION 2-5

6. A screen displays requesting additional information for this system. Select the **Model Type** of the system (See Table 2-1 for Model Type information) and enter the **System ID**. Leave the **Serial Number** field blank for all MR systems. See Illustration 2-6.

TABLE 2-1
MR MODEL TYPES

System	Software Version	Model Type(s)
MR450	20.1_MX_xxxx.x	MR
MR750	20.0_MX_xxxx.x	MR or Discovery_MR750
HDxt	15.0_MX_xxxx.x	MR or Signa_HDxt
HDi	15.0_MX_xxxx.x	MR or Signa_HDi
HDx	14.0_MX_xxxx.x	MR or SIGNA-HDx
HDe	HDE14.0_MX_xxxx.x	SIGNA-HDe
HD	12.0_MX_xxxx.x	MR or Excite12
Excite II	11.0_MX_xxxx.x	Excite11
Excite II – G3	G3.0_MX_xxxx.x	3TG3
Excite II – E2	E2.0_MX_xxxx.x	3TG3
Openspeed w/Excite	OPENMR4x_MX_xxxx.x	HFO4
Ovation w/Excite	OPENMR4_MX_xxxx.x	MFO4
Signa LX	9.1	MR or LX91
Goldseal	Various	MR

Add New System

The system you entered does not exist in our database. Enter the details for this system here to add it to our database, or close this window to re-enter the host ID

* Modality: MR

* Model Type: 

Serial Number:

* Host ID: 0384924

System ID:

USN:



ADD NEW SYSTEM POP-UP
ILLUSTRATION 2-6

7. The *Order Information* screen displays, with a new column for the added system. Items on the order that correspond to an option key for the selected Model Type have a checkbox next to it. See Illustration 2-7.

Order Information

FDO: 0002 002612 /
 Site: CHILDRENS HOSPITALS & CLINICS
 Billing Number: 000146914
 Address: ACCOUNTS PAYABLE
 Street: 2910 CENTRE POINT DR
 City, State Zip: ROSEVILLE MN 55113

Instructions
[Add A System](#)
[Generate Licenses](#)
[View System Information](#)

3. Generate Licenses key

Add a System **Generate Licenses**

Host : 0384924 Details	Description	Qty Ordered	Qty Used
System : 023456MR			
Modality : MR			
Model Type : Excite11			
☐ Check All			
	M3033KK	1	0
	M3090KJ	1	0
	S7502JS	1	0
☒	M3090PL	1	0
	M1090PM	1	0
	M1033PA	1	0
	M3033MV	1	0
	M80501DA	1	0
	S7501AF	1	0
	M3087N	1	0
	M1000NT	1	0
	M1087SP	1	0
	S7502NB	1	0
	M1000MN	1	0
	M1000MV	1	0

2. Check All button

1. Checkbox

GENERATING LICENSES FROM ORDER INFORMATION SCREEN
 ILLUSTRATION 2-7

8. Select all options that you want to generate for the system, or click **<Check All>** if there is only one system on the order, and then click **<Generate Licenses>**.
9. Read through the selected options, customer information, and system information. If everything is correct, click **<Confirm Selection>**. See Illustration 2-8.

Order Information

Site: CHILDRENS HOSPITALS & CLINICS
 FDO: 0002 002612 /
 Order Date: 2003-04-30
 Billing Number: 000146914
 Address: ACCOUNTS PAYABLE
 Street: 2910 CENTRE POINT DR
 City, State Zip: ROSEVILLE MN 55113

PRODUCT	ASSIGNED CATALOG ITEMS	
	Catalog No.	Description
Host : 0384924 System : 023456MR Modality : MR Model Type : Excite11	M3090PL	CONNECTPRO FOR EXCITE

Confirm Selection(s) **Back** **Confirm Selection button**

CONFIRM KEY GENERATION SCREEN
 ILLUSTRATION 2-8

10. The license information displays on the *Order Information* screen. A key displays in the System column of all generated options. Point the mouse to the key symbol to view the option string. See Illustration 2-9.

Order Information

FDO: 0002 002612 /
Site: CHILDRENS HOSPITALS & CLINICS
Billing Number: 000146914
Address: ACCOUNTS PAYABLE
Street: 2910 CENTRE POINT DR
City, State Zip: ROSEVILLE MN 55113

Instructions
[Add A System](#)
[Generate Licenses](#)
[View System Information](#)

Add a System

Generate Licenses

Host : 0384924 System : 023456MR Modality : MR Model Type : Excite11 Check All	Catalog Number	Description	Qty Ordered	Qty Used
			Quantity available & used	
	M3033KK	APA2 FOR FORWARD PRODUCTI	1	
	M3090KJ	2D FIESTA FOR EXCITE	1	0
	S7502JS	CONNECTPROPLUSBARCODE EX	1	0
	M3090PL	CONNECTPRO FOR EXCITE	1	1
	hisir permanent 0 1c7ba4953c13e6d470d00309	BAR CODE READER/CONNECTPRO	1	0
	M3033KK	PERFORMED PROCEDURE STEP	1	0
	M3033MVV	IDRIVE PRO PLUS FOR EXCIT	1	0
	M80501DA	OM-AW4.1-LINUX-USA	1	0
	S7501AF	DUAL ARRAY	1	0
	M3087N	NON OBSOLESCENCE	1	0

Option Generated Key


CONFIRM KEY GENERATION SCREEN ILLUSTRATION 2-9

2-1-3 Adding a System to an Order (For Upgrades)

- Go to the eLicense website <http://elicense.gehealthcare.com> and log on using your single sign-on username and password.
- From the eLicense homepage, enter the System ID, then select System ID and Modality of MR as shown in Illustration 2-10, then select **Go**. If the System ID is not associated with keys in eLicense, you can select Host ID from the pulldown and enter the Host ID of the original computer. To determine the Host ID of the system, open a C-shell and type **sysinfo**. Also, the original order number can be used if known (not the new order number).

Note:

It is very important that you locate all of your option keys prior to creating new keys on upgrades. Older system types like Signa LX do not have option keys in eLicense.

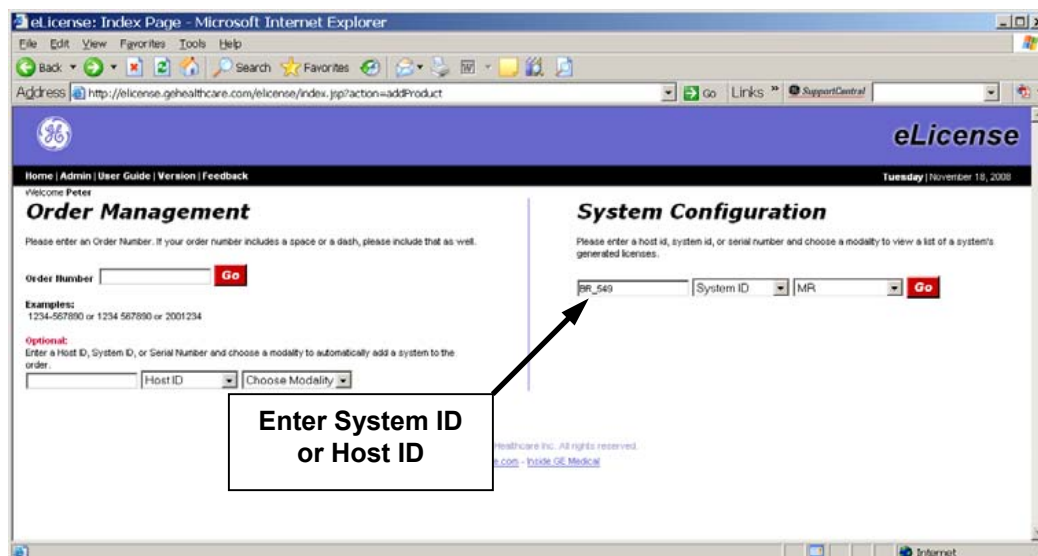

ORDER MANAGEMENT/SYSTEM CONFIGURATION SCREEN ILLUSTRATION 2-10

- Record the Order Number(s) as shown in Illustration 2-11.

Note

The Order ID column contains all of the order numbers of any key created for the system. It is good practice to record these numbers prior to doing the upgrade.

eLicense: System Information - Microsoft Internet Explorer

Address: <http://elicense.gehealthcare.com/elicense/search.do?action=searchResults>

System Information

Modality: MR
Model Type: Signa_HDxT
Serial Number:
Host ID: 973375730
System ID: KAISERM
USIT:

Instructions:
[Generate New Licenses](#)
[Email Licenses](#)
[Save Licenses](#)
[Update List](#)
[Reconfigure System](#)

Buttons: **Generate New Licenses** | **Email Licenses** | **Save Licenses** | **Print Licenses** | **Update List** | **Reconfigure**

Generate Date	Generate Type	Order ID	Genesis ID	Catalog ID	Description	Floatable	Application Name/Key/Number of Users	Generated By	Modified By	Modified Date	Download
2008-10-29	Generated From Order	2777625	2777625	M3033LZ	Fluoro Trigger MRA M3033LZ (14)	N	ftara permanent 0 3d06259b9b132393cf5a1361	212002418	212002418	2008-10-29	
2008-10-29	Generated From Order	2777625	2777625	M3090KJ	2D Fiesta M3090KJ (14)	N	fiesta2d permanent 0 3998249801137e6ae2cc0724	212002418	212002418	2008-10-29	
2008-10-29	Generated From Order	2777625	2777625	M3335KJ	Lava-XV M3335KJ (14)	N	lavaxv permanent 0 6c6985a0581342c233383061	212002418	212002418	2008-10-29	
2008-10-29	Generated From Order	2777625	2777625	M3033MV	iDrive Pro Plus M3033MV (14)	N	rtca permanent 0 300184ac6713499636f20082	212002418	212002418	2008-10-29	

RECORD THE ORDER NUMBER
ILLUSTRATION 2-11

- Select Edit System Information.

eLicense

Close Window

Enter System Information

Modality: MR

Model Type:

Serial Number:

Host ID:

MAC ID:

System ID:

USIT:

Continue

Annotations:
- Select Model Type (Table 2-1)
- Enter new Host ID if replacing host computer

ENTER SYSTEM INFORMATION
ILLUSTRATION 2-12

- Select the “Model Type” of the system (See Table 2-1 for Model Type Information).

Note


If the Host computer is being replaced, delete the Host ID in the Edit System Information GUI and change to the new one. Do not transfer the NIC card to the new computer. All previous option keys are automatically regenerated with this Host ID.

- Select **Continue** on the GUI
- Select **Generate Licenses** on the eLicense GUI

Note

Never select the **Reconfigure** button. It can erase all option keys if it is selected.

- Enter the “New / Upgrade Order Number” that was recorded in step 3.
- There should be additional Options noted as shown in Illustration 2-13.

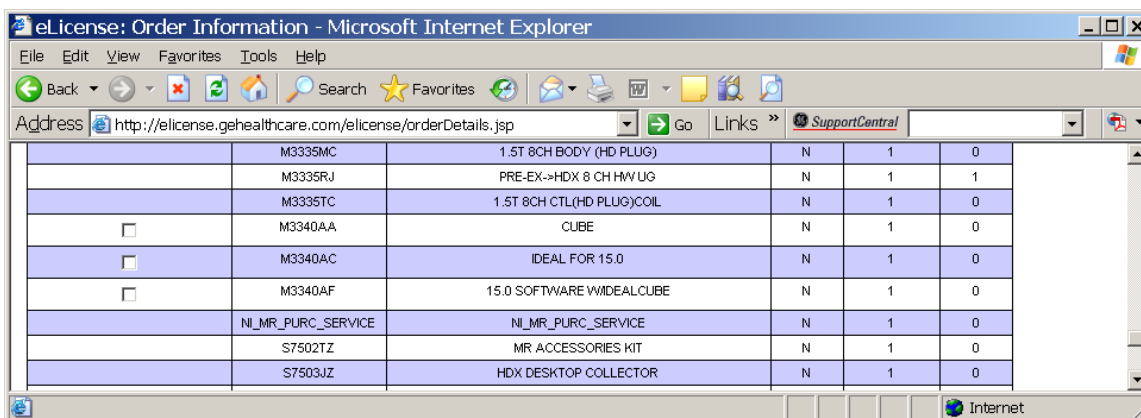

The screenshot shows a web browser window titled "eLicense: Order Information - Microsoft Internet Explorer". The address bar shows "http://elicense.gehealthcare.com/elicense/orderDetails.jsp". The main content area displays a table with the following data:

	M3335MC	1.5T 8CH BODY (HD PLUG)	N	1	0
	M3335RJ	PRE-EX->HDX 8 CH HW UG	N	1	1
	M3335TC	1.5T 8CH CTL(HD PLUG)COIL	N	1	0
☐	M3340AA	CUBE	N	1	0
☐	M3340AC	IDEAL FOR 15.0	N	1	0
☐	M3340AF	15.0 SOFTWARE WIDEALCUBE	N	1	0
	NI_MR_PURC_SERVICE	NI_MR_PURC_SERVICE	N	1	0
	S7502TZ	MR ACCESSORIES KIT	N	1	0
	S7503JZ	HDX DESKTOP COLLECTOR	N	1	0

ADDITIONAL OPTIONS
ILLUSTRATION 2-13

- Click on the new options, then select **Generate Licenses**.

2-2 Viewing existing licenses for a system

You can retrieve a list of licenses that have been previously generated for a system by searching from the eLicense home page.

- Go to the eLicense website <http://elicense.gehealthcare.com> and log on using your single sign-on username and password.
- From the eLicense home page, enter the System ID or Host ID and modality, in the *System Configuration* section and click **<Go>**. The *System Information* page lists all the options that have been generated for that system. Other information shown is *Generate Date*, *Order ID*, *License information*, and the SSO of the user who generated the option keys. See Illustration 2-14.

System Information

[Edit System Information](#)

Modality: MR
Model Type: Excite11
Serial Number:
Host ID: 0384924
System ID: 023456MR
USH:

Instructions

[Generate New Licenses](#)
[Email Licenses](#)
[Save Licenses](#)
[Update List](#)
[Reconfigure System](#)

Generate New Licenses		Email Licenses		Save Licenses		Print Licenses		Update List		Reconfigure	
Generate Date	Generate Type	OCF	CSS Order ID	Genesis ID	Catalog ID	Description	License Text			Generated By	Download
2003-10-07	Generated From Order	0002	002612		M3090PL	ConnecPro (Excite 11.0)	hisris permanent 0 1c7ba4953c13e6d470d00309			212012346	

SYSTEM INFORMATION SCREEN
ILLUSTRATION 2-14

- An alternative method for retrieving a list of existing licenses is to enter the *Global Order Number* in the *Order Management* section and click **<Go>**. As long as keys were created, you can find the System ID or Host ID that it is tied to. In the *Order Assignment* page, locate the proper System ID and select **Details**. See Illustration 2-15.

The screenshot shows the 'eLicense: Order Information - Microsoft Internet Explorer' window. The address bar shows 'http://elicense.gehealthcare.com/elicense/orderDetails.jsp'. The page title is 'eLicense'. The navigation bar includes 'Home | Admin | User Guide | Version | Feedback | Logout' and the date 'Monday | November 9, 2009'. The main section is titled 'Order Assignment'. It displays system information: FBO: 2782058, Header ID: 12219001, Site: CRISP REGIONAL HOSPITAL, Address: CRISP REGIONAL HOSPITAL, Street: 302 N SEVENTH ST, City, State, Zip: CORDELE, GA, 31015, Country: US. Below this is a table of licenses. The first row is highlighted, and a callout box points to the 'Details' button for that row. The table has columns: Host, System, Serial Number, Modality, Model Type, Catalog Number, Description, Floatable, Qty Ordered, and Qty Used. The first row shows Host: 976907602, System: 229276MR, Serial Number: , Modality: MR, Model Type: MR, Catalog Number: 2121199-12, Description: UIG KIT, NEW BLOWER BRKTS & CAP MNTS, Floatable: N, Qty Ordered: 1, Qty Used: 0.

ORDER ASSIGNMENT DETAILS SCREEN
ILLUSTRATION 2-15


2-3 Distributing Licenses

There are four methods for distributing license keys generated at the eLicense server: email, save, print, or download via inSite. You can access these options from the eLicense *System Information* page. Follow step 2-2 of this document to access the *System Information* page, or, click the **<Details>** button on the *Order Information* screen. See Illustration 2-16.

FDO: 8937 029477 / 00000017998
Site: GEN ELECT PORTUGUESA S A R L
Billing Number: 931999953
Address: AV DO FORTE, NO.4,FRACCAO 'F'
Street: 2795-197 CARNAXIDE
City, State Zip: PORTUGAL

[Add A System](#)
[Generate Licenses](#)
[View System Information](#)

Add a System		Generate Licenses			
Host : 76543210 Details System : Modality : NM Model Type : INTEGRA Check All		Catalog Number	Description	Qty Ordered	Qty Used
 ^cb+/wwpub5dfdwjz4r6+dx33x-gy-ax3A		Z80431X	VG 5/8 INT SYST W/INTEGRA	1	0
		H3000YB	MILLENNIUM VG 5/8 ACQ	1	0
		H3400LC	GENIE P&R CE KIT	1	0
		H2600WC	INTEGRA CSI OPTION	1	1
		E8505GA	RECT FLOOD PHANTOM 76-801	1	0
		H3602MC	BAR PHANTOM	1	0
		H3400ZF	ENT L ANGLU/AGE KIT-PORTUGES	1	0

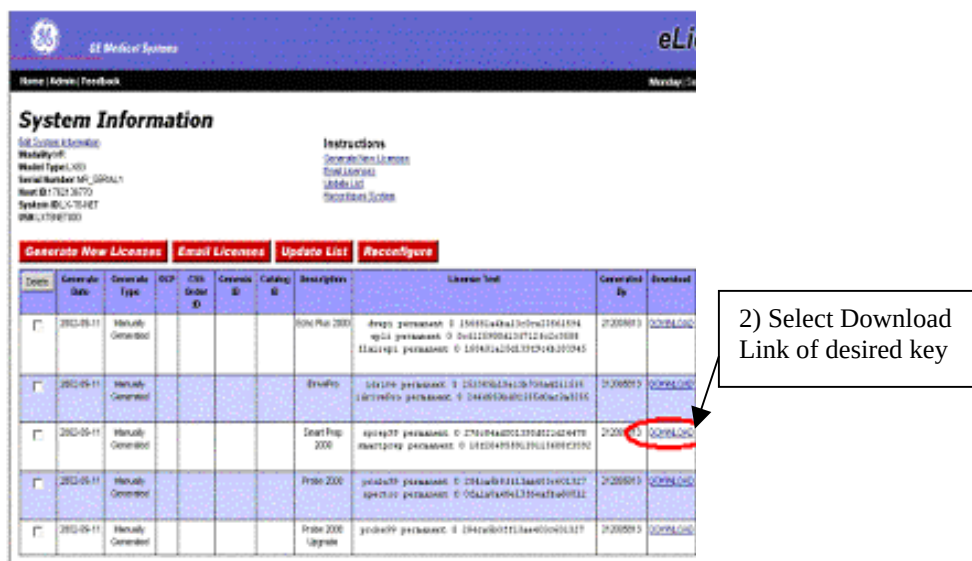

3. Wait for the request to be processed. A status page provides information on the downloading process.
4. Repeat steps 2-3 for all options that must be activated on the scanner.

2-3-2 Emailing Licenses

The license file can be emailed to any email address. Then, save this file to a floppy disk or CD, print it, or send it via FTP to the target system following the steps in Section 3.

1. From the *System Information* screen, select **[Email Licenses]**.
2. The *Email License Keys* window displays, as shown in Illustration 2-18.

EMAIL LICENSE KEYS POP-UP
ILLUSTRATION 2-18


3. Additional email can be added to this message.
 - a. If the recipient is a GE employee, enter their first and last name, and then click **<Search>**. Select the name you want from the list.
 - b. If the email address is an external address, type it in the **Address** field, and click **<Add>**.
4. To remove unwanted addresses, click **<Delete>**.
5. When all recipients are added, click **<Continue>**.
6. A verification window displays. Click **[Close Window]** in the upper right corner.

2-3-3 Saving Licenses to a Disk

You can save the license file to a floppy disk for installation via Guided Install on the scanner. The file can also be saved to the hard drive, and later transferred to a disk, printed, or sent via FTP to the scanner.

1. From the *System Information* screen, click **[Save Licenses]**.

2. On the *File Download* screen, select **[Save this file to disk]** and click **<OK>**. See Illustration 2-19.

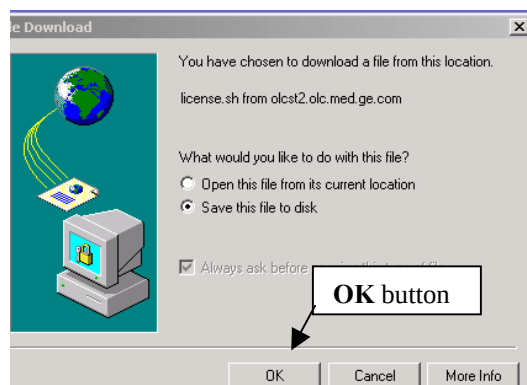


FILE DOWNLOAD SCREEN
ILLUSTRATION 2-19

3. At the **Save As** screen you can select the destination of the license file. See Illustration 2-20.
 - a. If saving to a disk, insert a disk into the computer's drive and select the drive from the menu.
 - b. If saving to the hard drive, select the folder to save the file to.
4. The file name must be **license.sh**. The software looks for this exact name.

Note

Windows will try to change the file name to license.txt. You must ensure that the file name is license.sh or installation will fail on the system.

5. In the *Save as Type* field, select **All Files**.
6. Click **<Save>** when finished.

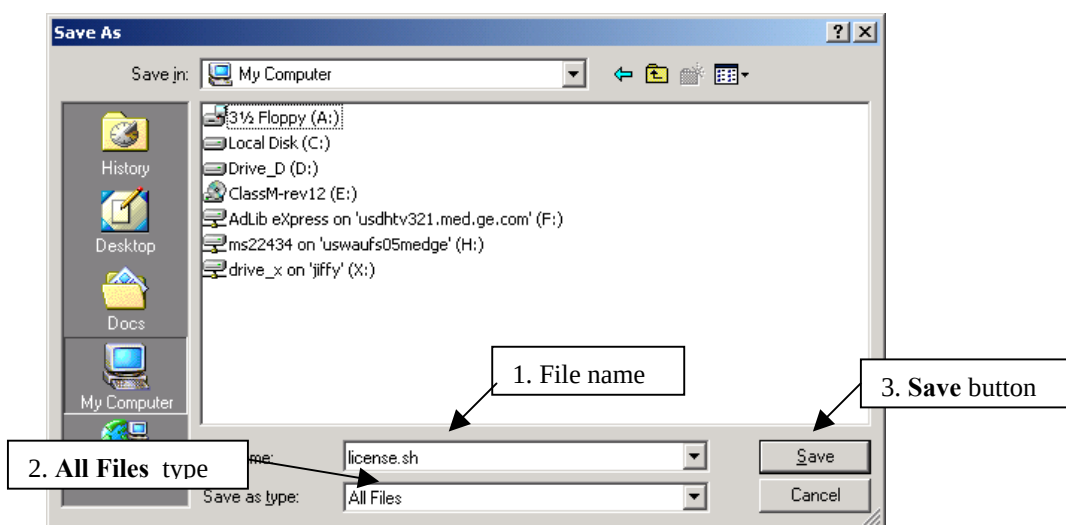

SAVE AS MENU
ILLUSTRATION 2-20

7. After the file is saved, a *Download Complete* window displays. Click **[Close]**.

System Information

[Edit System Information](#)

Modality: AWV
Model Type: AWG1
Serial Number:
Host ID: 80edb514
System ID: AWV_TEST_S1
USN:

Instructions
[Generate New Licenses](#)
[Email Licenses](#)
[Print Licenses](#)
[Update List](#)
[Reconfigure System](#)

Generate New Licenses			Email Licenses		Print Licenses		Update List		Reconfigure	
Delete	Generate Date	Generate Type	OCP	CSS Order ID	Genesis ID	Catalog ID	Description	License Text	Generated By	Download
☐	2002-09-25	Manual Sync				B79971NC	BMD SAW Only (Global)	BMD 2437625804	212022386	
☐	2002-09-25	Manual Sync				B79991SA	CardiacIQ (Global)	CardIQ 1585323284	212022386	
☐	2002-09-18	Manually Generated				B7997K	3D	3D 3853491896	305008005	

EDIT SYSTEM INFORMATION BUTTON
ILLUSTRATION 2-22

- A pop-up window displays to allow modification of the system information. See Illustration 2-23.

Note

All changes to system information and host ID are tracked based on your single sign-on ID. Any suspicious usage of this feature will be rigorously checked to verify the modification was due to a valid problem.

Enter System Information

Modality: MR

Model Type: EXCITE12

Serial Number:

Host ID:

System ID:

USN:

Continue

Copyright © 2002 GE Medical Systems Inc. All rights reserved.

EDIT SYSTEM INFORMATION SCREEN
ILLUSTRATION 2-23

Model Type	System Description
EXCITE11	all 1.5T and 3.0T 11x systems without G3 software
3TG3	only for early 11x systems with G3 software. All other 11x systems choose "EXCITE11"
EXCITE12 or MR	all 1.5T and 3.0T 12x (HD) systems
SIGNA-HDe	all 1.5T and 3.0T HDe systems
SIGNA-HDx or MR	all 1.5T and 3.0T HDx systems
Signa_HDxT or MR	all 1.5T and 3.0T HDxt systems
Signa_HDi or MR	all 1.5T Hdi systems
Discovery_MR750 or MR	MR750
MR	MR450

Note

For Upgrades: Change the **Model Type**, and if necessary, update the **Host ID** (see Section 2-1-1).

Note

If the Host Id has been changed, before and after clicking [Continue] verify that the *Application Name Key* number on the *System Information* page actually changed. Depending on how many keys are regenerated, it may take up to a minute before the page refreshes and you see new numbers. If the numbers do not update, you probably did not enter the *Host ID* correctly; go to step 2-4-2 and re-enter the information again.

3- OPTION INSTALLATION

There are five methods for installing eLicense option keys on the system. The details of each method are given in the following sections.

Note

The preferred method of installing option keys is to burn to a CD-ROM (see Section 3-2). Download may not install all option keys because dependent keys are not on the system.

IMPORTANT

WHEN INSTALLING OPTIONS ON A SYSTEM MAKE SURE THE OPTION MATCHES THE UPGRADE BEING PERFORMED. ALL OTHER CONFLICTING OPTION KEYS SHOULD BE REMOVED. FAILURE TO REMOVE THEM CAN RESULT IN A DE-RATED SYSTEM CONFIGURATION.

EXAMPLE 1: SYSTEM OPTIONS FOR HISPEED AND ECHOSPEED ARE BOTH INSTALLED. THE RESULT IS THE SYSTEM WILL PERFORM AS A HISPEED NOT AN ECHOSPEED.

TO RESOLVE – THE HISPEED OPTION MUST BE REMOVED TO ALLOW THE SYSTEM TO CONFIGURE AS AN ECHOSPEED.

EXAMPLE 2: SYSTEM OPTIONS FOR 8-CHANNEL AND 16-CHANNELS ARE BOTH INSTALLED. THE 8 CHANNEL OPTION MUST BE REMOVED TO ALLOW THE SYSTEM TO CONFIGURE AS A 16-CHANNEL.

3-1 eLicense Floppy

Note

On new systems the option floppy is no longer a supported feature. Verify the host PC has a floppy drive prior to creating an eLicense Floppy.

All of the eLicense option keys can be saved onto one floppy disk. After creating this floppy disk, follow these instructions to load them onto the system.


Note

The option key file must end with a **.sh** suffix. If the file was saved as anything other than **license.sh**, the file will not be recognized when loaded onto the scanner. Please rename the file to **license.sh**.

1. Insert the floppy disk into the host computer floppy drive.
2. Select the **[Service Tools Desktop]**.
3. Select **[Guided Install]** from the Service Desktop Manager and press **[Start...]** to open Guided Install.
4. A command window displays **"Type Root Password"**. Type the system's root password and press **<Enter>**. **"operator"** is the Signa default root password. If the customer has changed the root password, you need to obtain the correct password from the customer.
5. Once the Guided Install window opens, select the **Option Keys** tab. See Illustration 3-1.
6. Select the **eLicense Floppy** check box on the Install Option Key portion of the window, and then select **<Install>**.


Note

Select the **Floppy** checkbox when using an option floppy sent with the system and not generated through eLicense. The eLicense floppy will NOT work when the floppy checkbox is selected. Please use the **eLicense Floppy** selection for keys generated with eLicense, and **Floppy** for factory generated keys.

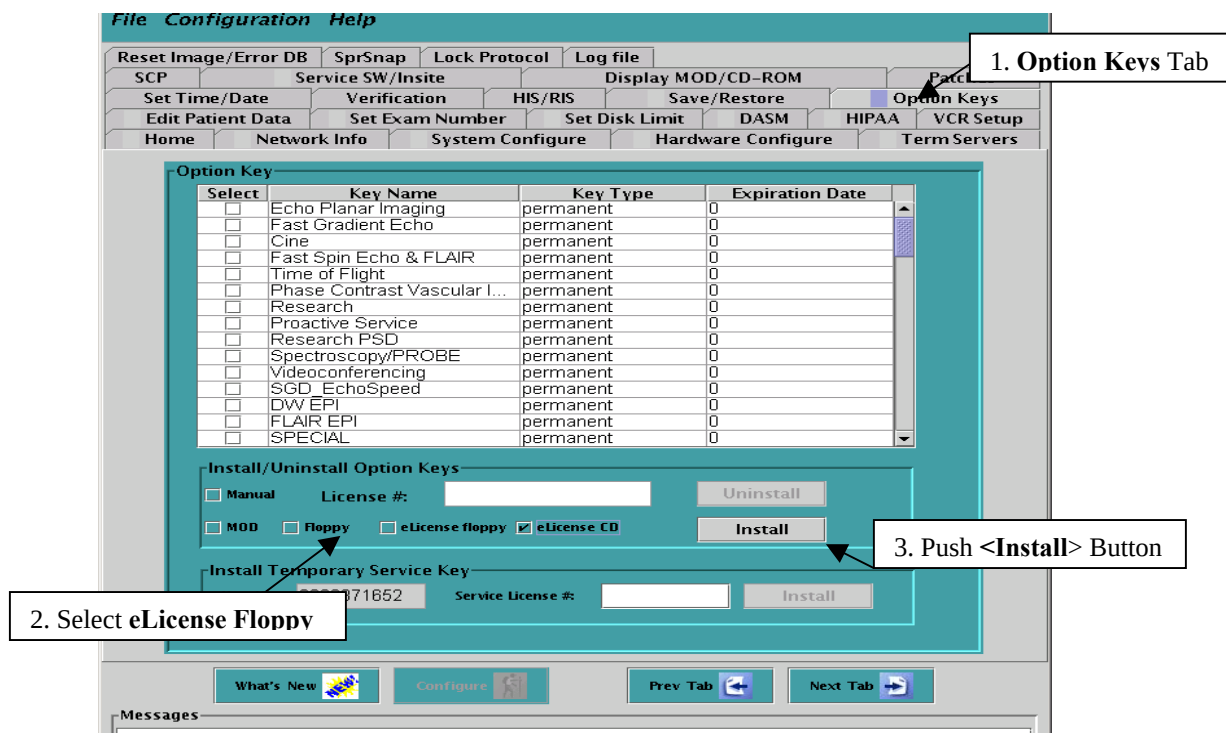

GUIDED INSTALL WINDOW
ILLUSTRATION 3-1

7. The option key file will begin loading. Load times per option key are very short. Typically less than one minute. When the keys are finished loading, a pop-up displays as shown in Illustration 3-2, listing all the installed options. Scantools and EchoPlus options are shown in the Illustration.

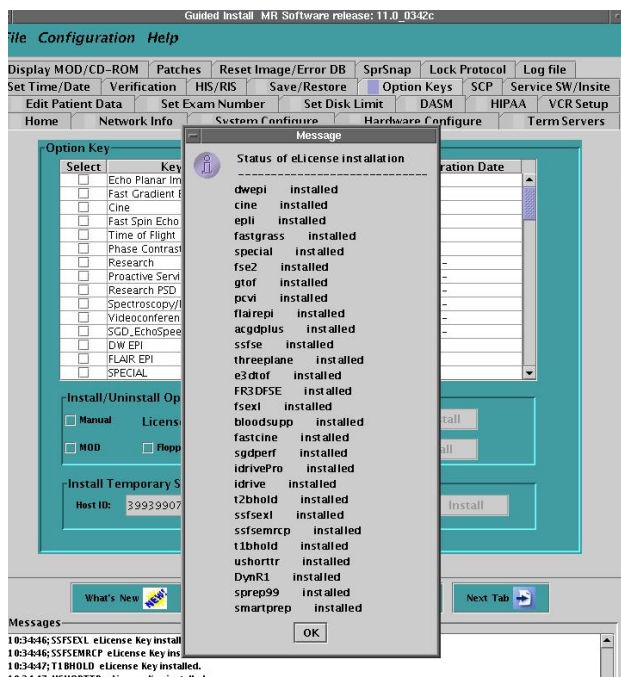

ELICENSE FLOPPY KEY INSTALLATION COMPLETE MESSAGE
ILLUSTRATION 3-2

8. Verify that the option list in the Guided Install GUI is updated with the options just loaded, if they were previously “Not Installed” or of type “FlexTrial”.

IMPORTANT

WHEN INSTALLING OPTIONS ON A SYSTEM MAKE SURE THE OPTION MATCHES THE UPGRADE BEING PERFORMED. ALL OTHER CONFLICTING OPTION KEYS SHOULD BE REMOVED. FAILURE TO REMOVE THEM CAN RESULT IN A DE-RATED SYSTEM CONFIGURATION.

EXAMPLE 1: SYSTEM OPTIONS FOR HISPEED AND ECHOSPEED ARE BOTH INSTALLED. THE RESULT IS THE SYSTEM WILL PERFORM AS A HISPEED NOT AN ECHOSPEED.

TO RESOLVE – THE HISPEED OPTION MUST BE REMOVED TO ALLOW THE SYSTEM TO CONFIGURE AS AN ECHOSPEED.

EXAMPLE 2: SYSTEM OPTIONS FOR 8-CHANNEL AND 16-CHANNELS ARE BOTH INSTALLED. THE 8 CHANNEL OPTION MUST BE REMOVED TO ALLOW THE SYSTEM TO CONFIGURE AS A 16-CHANNEL.

9. Select **File -> Quit** from the main menu bar to exit the Guided Install tool.
10. Reboot the computer, and verify that all options are available.
11. Complete a **Save Info** to retain all newly installed options on an Info MOD.

3-2 eLicense CD-ROM

All of the eLicense option keys can also be saved onto one CD-ROM. **This is the preferred method for installing option keys.** Once you create the CD, the options can be loaded onto the system, described in section 3-2-2.

3-2-1 Making the eLicense CD-ROM

Once the option key file is downloaded from the eLicense site to the hard drive of a computer, it can be burned onto a CD. The MR system requires that the disc be ISO-compliant. This format type cannot be burned with Windows XP's built-in disc burning. To burn an ISO-compliant disc, you must use a program such as Sonic RecordNow, which burns ISO-compliant discs by default.

1. Start Sonic RecordNow.

Note


To obtain a licensed copy of Sonic go to <ftp://globe/globepro/swstore/applications/HelpDesk/SONIC> and download the Sonic Record Now 7.exe file or call GE Healthcare's Help Desk.

2. From Sonic's home page, select the **Data Projects** tab, and then click **Data Disc** (see Illustration 3-2-1-1).



CREATING A DATA DISC
ILLUSTRATION 3-2-1-1

3. Click **Add Files and Folders** (see Illustration 3-2-1-2), and then browse your PC to add the **license.sh** to the **Data Disc** list.

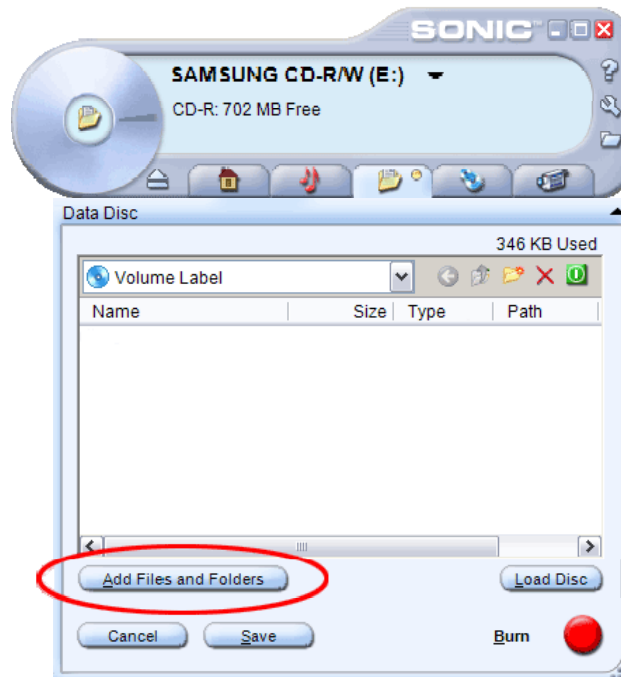


ADDING FILES
ILLUSTRATION 3-2-1-2

4. Put a blank CD in the computer's CD-ROM drive and click **<Burn>**.

3-2-2 Loading options from eCD

After creating this CD-ROM, follow these instructions to load the keys onto the system.

1. Insert the CD-ROM into the host computer CD drive.
2. Select the **[Service Tools Desktop]**.
3. Select **[Guided Install]** from the Service Desktop Manager and press **[Start...]** to open Guided Install.
4. A command window displays **"Type Root Password"**. Type the system's root password and press **<Enter>**. **"operator"** is the Signa default root password. If the customer has changed the root password, you need to obtain the correct password from the customer.
5. Once the Guided Install window opens, select the **Option Keys** tab. See Illustration 3-7.
6. Select the **eLicense CD** check box on the Install Option Key portion of the window, and then select **<Install>**.

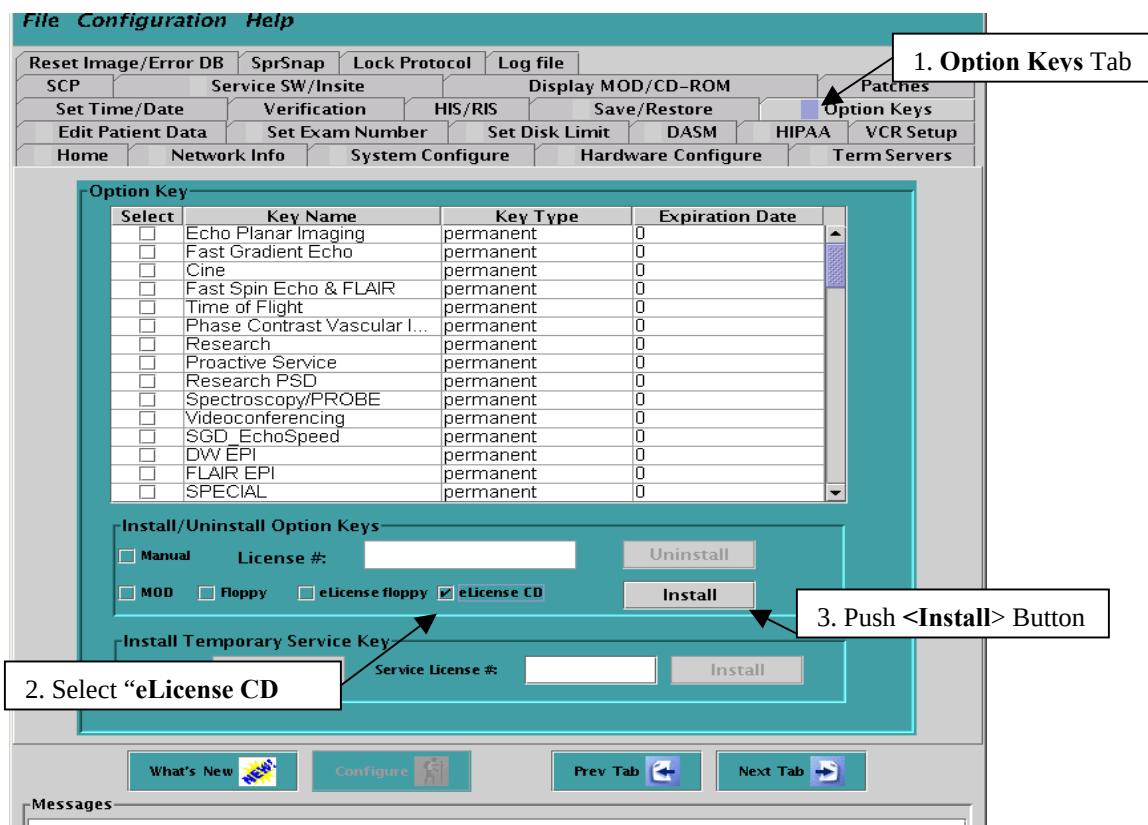


GUIDED INSTALL WINDOW
ILLUSTRATION 3-7

7. When the keys are finished loading, a pop-up displays the list of the installed options.
8. Verify the option list in the Guided Install GUI is updated with the options just loaded.
9. Select **File -> Quit** from the main menu bar to exit the Guided Install tool.
10. Reboot the computer, and verify that all options are available.
11. Complete a **Save Info** to retain all newly installed options on an Info MOD.

3-3 Option Key Auto Download

If the auto download procedure is performed as described in section 2-3-4, the keys will be on the scanner waiting for a reboot.

1. Reboot the MR scanner after all option keys have been downloaded from the eLicense server.
2. Verify that a pop-up displays the message **"New option keys have been installed"** after the **"Please do not interact with the system until this message disappears..."** window disappears. See Illustration 3-8 for an example of this pop-up.

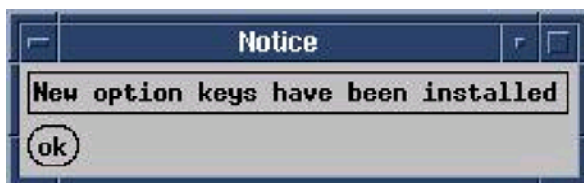

BOOT UP POP-UP MESSAGE AFTER INSTALLING NEW OPTIONS
ILLUSTRATION 3-8

- a. Open Guided Install, and verify the downloaded option keys are installed as permanent keys on the **Option Keys** tab.
3. Complete a **Save Info** to retain all newly installed options on an Info MOD.

3-4 FTP Option keys to scanner


If an inSite checkout has not been completed on the scanner, but the system is connected to the hospital network, the option keys can be sent to the scanner via FTP. The option key file must be saved to the computer hard drive before continuing.

1. Obtain the scanner IP address and connect to the system.
2. Send the option key file from the desktop to the scanner via FTP.
3. Move the file to the temp directory by typing "**mv license.sh /usr/tmp**".
4. Rename the file by typing "**mv license.sh flexkeys_install**".
5. Reboot the system.
6. Verify that a pop-up displays the message "**New option keys have been installed**" after the "**Do not interact with this system...**" message disappears.
 - a. Open Guided Install, and verify that the option keys are installed as permanent keys on the **Option Keys** tab.
7. Complete a **Save Info** to retain all newly installed options on an Info MOD.

3-5 Manual Option Key Entry

If all other methods of option key installation are unavailable, you can enter the keys manually. A printed copy of the keys is needed to transcribe the strings.

1. Select the **[Service Tools Desktop]**.
2. Select **[Install...]** from the Service Desktop Manager to open Guided Install.
3. A command window displays "**Type Root Password**". Type the system's root password and press **<Enter>**. The Signa default root password is "**operator**". If the customer has changed the root password, you need to obtain the correct password from the customer.
4. After the Guided Install window opens, select the **Option Keys** tab. See Illustration 3-9.
5. Select the **Manual** check box on the Install Option Key portion of the window, enter the 24-character license key string into the text box, and select **<Install>**. This field is case-sensitive.

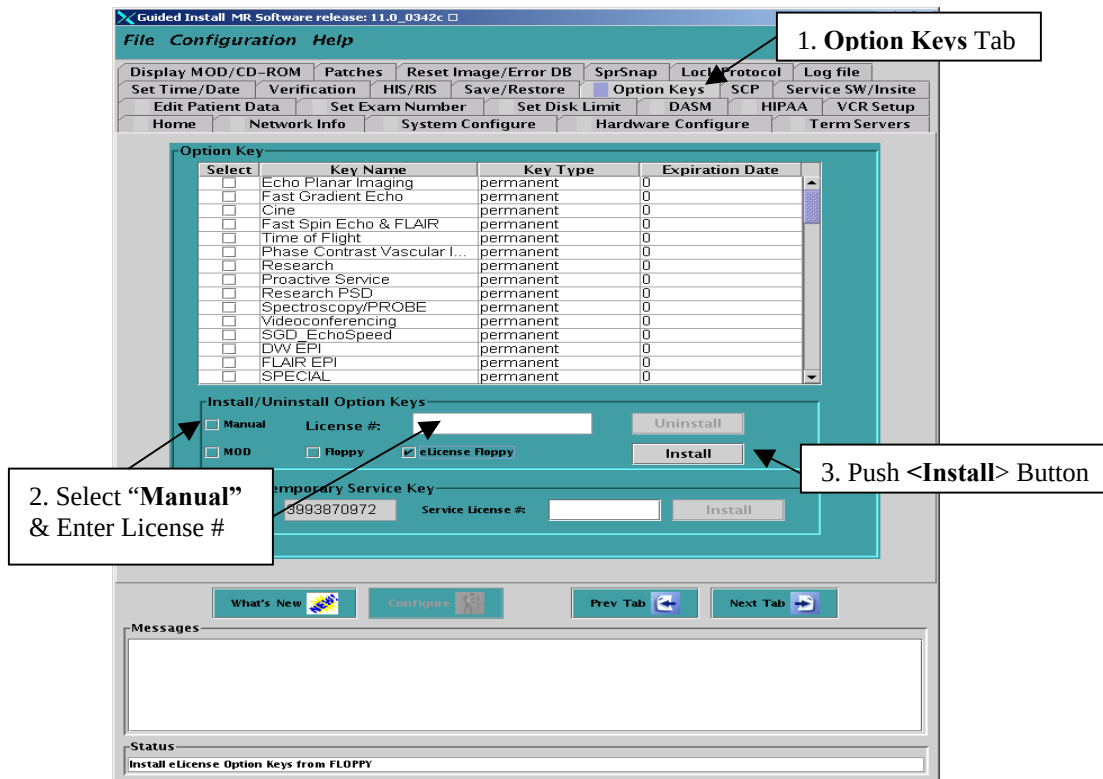

GUIDED INSTALL WINDOW
ILLUSTRATION 3-9

6. The option key will load immediately.
7. Verify the option list in the Guided Install GUI is updated with the option just loaded.
8. Repeat steps 5-7 for all the option keys that you are installing.
9. Select **File -> Quit** from the main menu bar to exit the Guided Install tool.
10. Reboot the computer, and verify that all options are available.
11. Complete a **Save Info** to retain all newly installed options on an Info MOD.

4- OPTION KEY CATALOG MATRIX

Refer to Table 4-1 for a list of available electronic options for MR Scanners.

Note

The Options Activated list below may change periodically. To see the most current keys, refer to the official drawing prints for each part.

TABLE 4-1
OPTION KEY CATALOG MATRIX

Option Catalog	Options Activated	Customer Protocol
L1000EE EXCITE11 GS CLASS 1.5T ES+EX 4CHN L1000EF EXCITE11 GS CLASS 1.5TES+EX 8CHF L1000EE EXCITE11 GSC MOB 1.5T ES+EX 4 CHN L1000EF EXCITE11 GS CLASS 1.5TES+EX 8CHM L2000EE EXCITE11 GOLDSEAL 1.5T EXCITE 4CH L2000EF EXCITE11 GOLDSEAL 1.5T EXCITE 8CH L3000EE EXCITE11 GOLDSEAL 1.5T EXCITE 4CH L3000EF EXCITE11 GOLDSEAL 1.5T EXCITE 8CH	▪ acgdplus ▪ asset ▪ asset91 ▪ bloodsupp ▪ cine ▪ dwepi ▪ DynR1 ▪ e3dtof ▪ epli ▪ fastcine ▪ fastgrass ▪ fiesta2d ▪ fiesta3d ▪ fiesta3dfs ▪ fiestac ▪ flairepi ▪ FR3DFSE ▪ fse2 ▪ fsexl ▪ ftmra ▪ gtof ▪ idrive ▪ idrivePro ▪ pcvi ▪ sgd ▪ smartprep ▪ special ▪ sprep99 ▪ ssfse ▪ ssfsemrcp ▪ ssfsexl ▪ t1bhold ▪ t2bhold ▪ threeplane ▪ ushorttr	▪ ACGD Plus ▪ Asset ▪ Asset Plus ▪ Blood Suppression ▪ Cine ▪ Diffusion Weighted epi ▪ Dynamic R1 ▪ 3D Time of Flight ▪ Echo Planer Imaging ▪ Fast Cine ▪ Fast Gradient Echo ▪ 2D Fiesta ▪ 3D Fiesta ▪ 3D Fast SAT Fiesta ▪ Phased Cycle Fiesta ▪ Flair EPI ▪ 3D FSE ▪ Fast Spin Echo & Flair ▪ FSE-XL ▪ Fluoro Triggered MRA ▪ Time of Flight ▪ iDrive ▪ iDrive Pro ▪ Phase Contrast Vascular Imaging ▪ SGD_EchoSpeed ▪ Smart Prep ▪ SPECIAL ▪ Smart Prep 2000 upgrade ▪ SSFSE ▪ SSFSE MRCP/Body Pak ▪ MRCP3 ▪ T1 Breath Hold ▪ T2 Breath Hold ▪ Three Plane Localizer ▪ Ultra Short TR

Option Catalog	Options Activated	Customer Protocol
M3000LL SIGNA-HDx E-III 4CH To E-III 8CH Upgrade M3335PB SIGNA-HDx 1.5T LCC Magnet HDx M3335PD SIGNA-HDx 1.5T LCC Twin Magnet HDx M3335PF SIGNA-HDx 3.0T G3 Magnet HDx	channel8	8-Channel
M3000PE SIGNA-HDx 16 Channel M3335JA SIGNA-HDx 16 Channel HDx 1.5T System M3335JG SIGNA-HDx 16 Channel Twinspeed HDx 1.5T System M3335JK SIGNA-HDx 16 Channel Twinspeed HDx 3T System	channel16	16-Channel
M3033LJ SIGNA-HDx INTEGRATED PATIENT COMFORT MODULE	ipcm	Integrated Patient Comfort Module
M3033MS ScanTools 11.0	- cine - epli - fastgrass - fse2 - pcvi - gtof - acgdplus - flairepi - special - ssfse - threeplane - e3dtof - fsexl - bloodsupp - fastcine - idrive - t2bhold - ssfsemrcp - t1bhold - ushorttr - ssfsexl - DynR1 - smartprep	- Cine - Echo Planer Imaging - Fast Gradient Echo - Fast Spin Echo & Flair - Phase Contrast Vascular Imaging - Time of Flight - ACGD Plus - Flair EPI - SPECIAL - SSFSE - Three Plane Localizer 3D Time of Flight - FSE-XL - Blood Suppression - Fast Cine - iDrive - T2 Breath Hold - SSFSE MRCP/Body Pak - T1 Breath Hold - Ultra Short TR - MRCP3 - Dynamic R1 - Smart Prep

Option Catalog	Options Activated	Customer Protocol
M3033KS ASSET for EXCITE 1.5T	- Asset - asset91	- ASSET - Asset Plus
M3033NS EXCITE12 Vibrant Bilateral Breast	Breast2	VIBRANT
M3335LT SIGNA-HDx Vibrant Bilateral Breast		
M3033NJ PROSE for EXCITE	Prose	Probe 3D Prostate
M3033MB FIESTA-C	Fiestac	Phased Cycled Fiesta
M3033LW 3D FRFSE for EXCITE	FR3DFSE	3D FSE
M3033MT iDrive Pro for Excite	idrivePro	iDrive Pro
M3033MW iDrive Pro Plus 1.5T	rtca	iDrive Pro Plus
M3090LP EchoPlus	dwepi	Diffusion Weighted epi
M3335KM HDx MR Echo	- realcard - tmcourse - delenhmt	- Fiesta Wall-Motion - Timecourse - MDE
M3033PJ Performed Proc Step for EXCITE	Pps	PPS
M3090PL ConnectPro for EXCITE	Hisris	Modality WorkList
M3033MJ VIBRANT Bilateral Breast	breast2	VIBRANT
M3033ND VIBRANT-II	breast2	VIBRANT-II
M3335KK BREASE	breastsvq	BREASE
M3033NC LAVA	Lava	Lava
M3335JS SIGNA-HDx Multi-Phase (variable delays)	MphVar	Multi-phase Variable Sequences
M3335KJ LAVA-XV	lavaxv	LAVA-XV
M3335KJ SIGNA-HDx Lava-XV		
M3033MA PROPELLER	- propt2 - propdwi	Propeller

Option Catalog	Options Activated	Customer Protocol
M3100LA SIGNA-HDx 1.5T MNS BASE KIT M3100LD SIGNA-HDx 1.5T MNS UPGRADE KIT M3334WE SIGNA-HDx 3.0T 8KW MNS BASE OPTION KIT M3100WH SIGNA-HDx E-III 3T MNS W-O CON CAB PRESENT	▪ mns ▪ sage ▪ probep ▪ probe3d ▪ spectro ▪ probe99	▪ Multi-Nuclear Spect ▪ Spectro Analysis GE ▪ Probe-P ▪ 3D Probe ▪ Spectroscopy/PROBE ▪ Probe 2000 Upg
M3100LE SIGNA-HDx Mutiple MRS M3334WJ SIGNA-HDx Mutiple MRS 3T	▪ probep ▪ probe3d ▪ spectro ▪ probe99	▪ Probe-P ▪ 3D Probe ▪ Spectroscopy/PROBE ▪ Probe 2000 Upg
M3033NH PROPELLER Brain Suite	▪ propdwi ▪ propt2	▪ Propeller DWI ▪ Propeller T2 ▪ Propeller T2 Flair
M3335KY T2 Mapping	▪ T2map	▪ T2 Mapping - Cartilage Integrity Package
M1033BT Brainwave Real Time	▪ brnwavrt	▪ Brainwave RT
M1033BS Brainwave Post Acquisition	▪ brnwavpa	▪ Brainwave PA
M1033BL BrainWave Light M3335MJ SIGNA-HDx BrainWave HW for HDx	▪ brnwavhw	▪ Brainwave Light Hardware
M3335KR BRAINWAVE 2.0 RT	▪ brnwavrt	▪ BrainWave 2.0 RT
M3335KS BRAINWAVE 2.0 PA - Console	▪ brnwavpa	▪ BrainWave 2.0 PA - Console
M3335KT BRAINWAVE FUSION - Console	▪ brnwavrt	▪ BrainWave Fusion - Console
M3090KK 3D FIESTA for EXCITE	▪ fiesta3d	▪ 3D FIESTA
M3335KW Cosmic 3D M3335KW SIGNA-HDx 3D Cosmic	▪ cosmic3d	▪ Cosmic 3D

Option Catalog	Options Activated	Customer Protocol
M3335PL SIGNA-HDx 1.5T Excite to HD to HDx cabinet UG components M3335PP SIGNA-HDx 1.5T HD to HDx Cabinet UG components	▪ acgdplus ▪ bloodsupp ▪ cine ▪ DynR1 ▪ e3dtof ▪ epli ▪ fastcine ▪ fastgrass ▪ flairpi ▪ FR3DFSE ▪ fse2 ▪ fsexl ▪ gtof ▪ idrive ▪ idrivePro ▪ pcvi ▪ sgd ▪ sgdperf ▪ smartprep ▪ special ▪ spre99 ▪ ssfse ▪ ssfsemrcp ▪ ssfsexl ▪ t1bhold ▪ t2bhold ▪ threeplane ▪ ushorttr ▪ lava ▪ dwepi ▪ asset ▪ asset91 ▪ merge2d ▪ bravo ▪ channel8 ▪ tmcourse ▪ delenhmt	▪ ACGD Plus ▪ Blood Suppression ▪ Cine ▪ Dynamic R1 ▪ 3D Time of Flight ▪ Echo Planer Imaging ▪ Fast Cine ▪ Fast Gradient Echo ▪ Flair EPI ▪ 3D FSE ▪ Fast Spin Echo & Flair ▪ FSE-XL ▪ Time of Flight ▪ iDrive ▪ iDrive Pro ▪ Phase Contrast Vascular Imaging ▪ SGD_EchoSpeed ▪ SGDPerf ▪ Smart Prep ▪ SPECIAL ▪ Smart Prep 2000 upgrade ▪ SSFSE ▪ SSFSE MRCP/Body Pak ▪ MRCP3 ▪ T1 Breath Hold ▪ T2 Breath Hold ▪ Three Plane Localizer ▪ Ultra Short TR ▪ Lava ▪ Diffusion Weighted epi ▪ Asset ▪ Asset Plus ▪ 2D Merge ▪ Bravo ▪ 8-Channel ▪ Time Course ▪ Delayed Enhancement

Option Catalog	Options Activated	Customer Protocol
M3335PM SIGNA-HDx 3.0T Excite to HD to HDx cabinet UG components M3335RW SIGNA-HDx G3 to HDX UG Fixed M3335RY SIGNA-HDx G3 to HDX UG T-R	▪ acgdplus ▪ bloodsupp ▪ cine ▪ e3dtof ▪ epli ▪ fastcine ▪ fastgrass ▪ flairrepi ▪ FR3DFSE ▪ fse2 ▪ fsexl ▪ gtof ▪ idrive ▪ idrivePro ▪ pcvi ▪ sgd ▪ sgdperf ▪ smartprep ▪ special ▪ sprepp99 ▪ ssfse ▪ ssfsemrcp ▪ ssfsexl ▪ t1bhold ▪ t2bhold ▪ threeplane ▪ ushorttr ▪ lava ▪ dwepi ▪ asset ▪ asset91 ▪ merge2d ▪ bravo ▪ channel8	▪ ACGD Plus ▪ Blood Suppression ▪ Cine ▪ 3D Time of Flight ▪ Echo Planer Imaging ▪ Fast Cine ▪ Fast Gradient Echo ▪ Flair EPI ▪ 3D FSE ▪ Fast Spin Echo & Flair ▪ FSE-XL ▪ Time of Flight ▪ iDrive ▪ iDrive Pro ▪ Phase Contrast Vascular Imaging ▪ SGD_EchoSpeed ▪ SGDPerf ▪ Smart Prep ▪ SPECIAL ▪ Smart Prep 2000 upgrade ▪ SSFSE ▪ SSFSE MRCP/Body Pak ▪ MRCP3 ▪ T1 Breath Hold ▪ T2 Breath Hold ▪ Three Plane Localizer ▪ Ultra Short TR ▪ Lava ▪ Diffusion Weighted epi ▪ Asset ▪ Asset Plus ▪ 2D Merge ▪ Bravo ▪ 8-Channel

Option Catalog	Options Activated	Customer Protocol
M3335PS SIGNA-HDx 3.0T HD to HDx cabinet UG components	▪ acgdplus ▪ bloodsupp ▪ cine ▪ e3dtof ▪ epli ▪ fastcine ▪ fastgrass ▪ flairapi ▪ FR3DFSE ▪ fse2 ▪ fsexl ▪ gtof ▪ idrive ▪ idrivePro ▪ pcvi ▪ sgd ▪ sgdperf ▪ smartprep ▪ special ▪ sprep99 ▪ ssfse ▪ ssfsemrcp ▪ ssfsexl ▪ t1bhold ▪ t2bhold ▪ threeplane ▪ ushorttr ▪ lava ▪ dwepi ▪ asset ▪ asset91 ▪ merge2d ▪ bravo	▪ ACGD Plus ▪ Blood Suppression ▪ Cine ▪ 3D Time of Flight ▪ Echo Planer Imaging ▪ Fast Cine ▪ Fast Gradient Echo ▪ Flair EPI ▪ 3D FSE ▪ Fast Spin Echo & Flair ▪ FSE-XL ▪ Time of Flight ▪ iDrive ▪ iDrive Pro ▪ Phase Contrast Vascular Imaging ▪ SGD_EchoSpeed ▪ SGDPerf ▪ Smart Prep ▪ SPECIAL ▪ Smart Prep 2000 upgrade ▪ SSFSE ▪ SSFSE MRCP/Body Pak ▪ MRCP3 ▪ T1 Breath Hold ▪ T2 Breath Hold ▪ Three Plane Localizer ▪ Ultra Short TR ▪ Lava ▪ Diffusion Weighted epi ▪ Asset ▪ Asset Plus ▪ 2D Merge ▪ Bravo

Option Catalog	Options Activated	Customer Protocol
M3335RD SIGNA-HDx EX to HDX 1.5 8CH HW UG F		
M3335RE SIGNA-HDx EX to HDX 1.5 8CH HW UG M	▪ acgdplus ▪ bloodsupp ▪ cine ▪ DynR1 ▪ e3dtof ▪ epli ▪ fastcine ▪ fastgrass ▪ flairepi ▪ FR3DFSE ▪ fse2 ▪ fsexl ▪ gtof ▪ idrive ▪ idrivePro ▪ pcvi ▪ sgd ▪ sgdperf ▪ smartprep ▪ special ▪ spre99 ▪ ssfse ▪ ssfsemrcp ▪ ssfsexl ▪ t1bhold ▪ t2bhold ▪ threeplane ▪ ushorttr ▪ lava ▪ dwepi ▪ asset ▪ asset91 ▪ merge2d ▪ bravo ▪ channel8	▪ ACGD Plus ▪ Blood Suppression ▪ Cine ▪ Dynamic R1 ▪ 3D Time of Flight ▪ Echo Planer Imaging ▪ Fast Cine ▪ Fast Gradient Echo ▪ Flair EPI ▪ 3D FSE ▪ Fast Spin Echo & Flair ▪ FSE-XL ▪ Time of Flight ▪ iDrive ▪ iDrive Pro ▪ Phase Contrast Vascular Imaging ▪ SGD_EchoSpeed ▪ SGDPperf ▪ Smart Prep ▪ SPECIAL ▪ Smart Prep 2000 upgrade ▪ SSFSE ▪ SSFSE MRCP/Body Pak ▪ MRCP3 ▪ T1 Breath Hold ▪ T2 Breath Hold ▪ Three Plane Localizer ▪ Ultra Short TR ▪ Lava ▪ Diffusion Weighted epi ▪ Asset ▪ Asset Plus ▪ 2D Merge ▪ Bravo ▪ 8-Channel
M3335RF SIGNA-HDx EX to HDX 1.5 8CH HW T-R		
M3335RG SIGNA-HDx EX to HDX 1.5 NONLCC 8CH UG		
M3335RJ SIGNA-HDx PRE-EX to HDX 8 CH HW UG		
M3335RL SIGNA-HDx CVI PRE-EX to HDX 8 CH UG		
M3335RM SIGNA-HDx TS PRE-EX to HDX 8 CH UG		
M3335RL SIGNA-HDx CVI PRE-EX to HDX 8 CH UG		
M3335RP SIGNA-HDx LX-5.5 to HDX UG F		
M3335RR SIGNA-HDx LX-5.5 to HDX UG T-R		
M3335RS SIGNA-HDx 4x-5x to HDX UG F		
M3335RT SIGNA-HDx 4x-5x to HDX UG T- R		

Option Catalog	Options Activated	Customer Protocol
M3335RK SIGNA-HDx PRE-EX to HDX 16 CH HW UG M3335RN SIGNA-HDx TS PRE-EX to HDX 16CH UG	▪ acgdplus ▪ bloodsupp ▪ cine ▪ DynR1 ▪ e3dtof ▪ epli ▪ fastcine ▪ fastgrass ▪ flairapi ▪ FR3DFSE ▪ fse2 ▪ fsexl ▪ gtof ▪ idrive ▪ idrivePro ▪ pcvi ▪ sgd ▪ sgdperf ▪ smartprep ▪ special ▪ sprep99 ▪ ssfse ▪ ssfsemrcp ▪ ssfsexl ▪ t1bhold ▪ t2bhold ▪ threeplane ▪ ushorttr ▪ lava ▪ dwepi ▪ asset ▪ asset91 ▪ merge2d ▪ bravo ▪ channel16	▪ ACGD Plus ▪ Blood Suppression ▪ Cine ▪ Dynamic R1 ▪ 3D Time of Flight ▪ Echo Planer Imaging ▪ Fast Cine ▪ Fast Gradient Echo ▪ Flair EPI ▪ 3D FSE ▪ Fast Spin Echo & Flair ▪ FSE-XL ▪ Time of Flight ▪ iDrive ▪ iDrive Pro ▪ Phase Contrast Vascular Imaging ▪ SGD_EchoSpeed ▪ SGDPerf ▪ Smart Prep ▪ SPECIAL ▪ Smart Prep 2000 upgrade ▪ SSFSE ▪ SSFSE MRCP/Body Pak ▪ MRCP3 ▪ T1 Breath Hold ▪ T2 Breath Hold ▪ Three Plane Localizer ▪ Ultra Short TR ▪ Lava ▪ Diffusion Weighted epi ▪ Asset ▪ Asset Plus ▪ 2D Merge ▪ Bravo ▪ 16-Channel

Option Catalog	Options Activated	Customer Protocol
M3335RZ SIGNA-HDx 3T94 E2 to HDX 16 CH UG M3335S SIGNA-HDx 3T94 E2 to HDX 16 CH UG	▪ acgdplus ▪ bloodsupp ▪ cine ▪ e3dtof ▪ epli ▪ fastcine ▪ fastgrass ▪ flairrepi ▪ FR3DFSE ▪ fse2 ▪ fsexl ▪ gtof ▪ idrive ▪ idrivePro ▪ pcvi ▪ sgd ▪ sgdperf ▪ smartprep ▪ special ▪ spre99 ▪ ssfse ▪ ssfsemrcp ▪ ssfsexl ▪ t1bhold ▪ t2bhold ▪ threeplane ▪ ushorttr ▪ lava ▪ dwepi ▪ asset ▪ asset91 ▪ merge2d ▪ bravo ▪ channel16	▪ ACGD Plus ▪ Blood Suppression ▪ Cine ▪ 3D Time of Flight ▪ Echo Planer Imaging ▪ Fast Cine ▪ Fast Gradient Echo ▪ Flair EPI ▪ 3D FSE ▪ Fast Spin Echo & Flair ▪ FSE-XL ▪ Time of Flight ▪ iDrive ▪ iDrive Pro ▪ Phase Contrast Vascular Imaging ▪ SGD_EchoSpeed ▪ SGDPerf ▪ Smart Prep ▪ SPECIAL ▪ Smart Prep 2000 upgrade ▪ SSFSE ▪ SSFSE MRCP/Body Pak ▪ MRCP3 ▪ T1 Breath Hold ▪ T2 Breath Hold ▪ Three Plane Localizer ▪ Ultra Short TR ▪ Lava ▪ Diffusion Weighted epi ▪ Asset ▪ Asset Plus ▪ 2D Merge ▪ Bravo ▪ 16-Channel
M3090LR Diffusion Tensor EXCITE	▪ tensor	▪ Diffusion Tensor
M3033NK Tractography	▪ tracto	▪ Tractography (FiberTrak)
M3090MA Probe 2000	▪ probe99 ▪ probep	▪ Spectroscopy/PROBE ▪ ProbePRESS
M3090MB Probe 2000 Upg for EXCITE	▪ probe99	▪ Probe 2000 Upg
M3090MC Sage 7 for EXCITE	▪ Sage	▪ Spectro Analysis GE
M3090MD Probe 3D Brain for EXCITE	▪ probe3d	▪ 3D Probe
M3033ME Navigator 3D Cardiac EXI	▪ navigator	▪ Navigator
M3033MD 3D Delayed Enhancement	▪ delenhmt3d	▪ 3D Delayed Enhancement
M3033MF 3D FatSat FIESTA Cardiac	▪ fiesta3dfs	▪ 3D FatSat FIESTA
M3033MH TRICKS	▪ tricks	▪ Tricks
M3090KJ 2D FIESTA for EXCITE	▪ fiesta2d	▪ 2D Fiesta

Option Catalog	Options Activated	Customer Protocol
M3033NG RealTime FIESTA	▪ realcard	▪ RealCARD
M3033NL 2D FatSat FIESTA	▪ fiesta2dfs	▪ 2D FatSat FIESTA
M3033NR 3D Flair	▪ flair3d	▪ 3D Flair
M3033MP FGRE-ET	▪ perfusion	▪ Perfusion
M3033MR Real-Time FGRE-ET	▪ Rtet	▪ Real-Time FGRE-ET
M3090KR 2D Delayed Enhancement EX	▪ delenhmt	▪ Delayed Enhancement
M3033MM Spiral Hi-Res for EXCITE	▪ spiralHR	▪ Spiral Hi-Res
M3033MZ Real-Time Spiral EXCITE	▪ spiralRT	▪ Spiral Real-Time
M3033MY Cardiac Tagging for EXCITE	▪ tagging	▪ Cardiac Tagging
M3033LY SmartPrep 2000 for EXCITE	▪ spre99	▪ Smart Prep 2000 upgrade
M3033LZ Fluoro Triggered MRA EXCITE	▪ ftmra	▪ Fluoro Triggered MRA
M3335KP MRA SWIFT	▪ swift	▪ MRA SWIFT
M3090LF Hightspeed Plus (Excite)	▪ hispeed	▪ High Speed
M3090LG Echospeed Plus (Excite)	▪ Sgd	▪ SGD_EchoSpeed
M3090LA Excite II MNS Kit	▪ mns ▪ spectro ▪ probe99 ▪ sage	▪ Multi-Nuclear Spect ▪ Spectroscopy/PROBE ▪ Probe 2000 Upg ▪ Spectro Analysis GE
M3090LB Excite II Fwd Prod MNS Kit	▪ mns ▪ spectro ▪ probe99 ▪ sage	▪ Multi-Nuclear Spect ▪ Spectroscopy/PROBE ▪ Probe 2000 Upg ▪ Spectro Analysis GE
M3090LC Excite II Legacy MNS upgrade	▪ mns ▪ spectro ▪ probe99 ▪ sage	▪ Multi-Nuclear Spect ▪ Spectroscopy/PROBE ▪ Probe 2000 Upg ▪ Spectro Analysis GE

Option Catalog	Options Activated	Customer Protocol
M3033NA 1.5T ScanTools 12.0	▪ cine ▪ epli ▪ fastgrass ▪ fse2 ▪ pcvi ▪ gtof ▪ acgdplus ▪ flairapi ▪ special ▪ ssfse ▪ threeplane ▪ e3dtof ▪ fsexl ▪ bloodsupp ▪ fastcine ▪ idrive ▪ t2bhold ▪ ssfsemrcp ▪ t1bhold ▪ ushorttr ▪ ssfsexl ▪ DynR1 ▪ smartprep ▪ sgd ▪ sgdperf ▪ idrivePro ▪ FR3DFSE ▪ sprep99	▪ Cine ▪ Echo Planar Imaging ▪ Fast Gradient Echo ▪ Fast Spin Echo & FLAIR ▪ Phase Contrast Vascular Imaging ▪ Time of Flight ▪ ACGD PLUS ▪ FLAIR EPI ▪ SPECIAL ▪ SSFSE ▪ Three Plane Localizer ▪ 3D Time of Flight ▪ FSE-XL ▪ Blood Suppression ▪ Fast Cine ▪ iDrive ▪ T2 Breath Hold ▪ SSFSE MRCP/Body Pak ▪ T1 Breath Hold ▪ Ultrashort TR/TE ▪ MRCP3 ▪ Dynamic R1 ▪ Smart Prep ▪ SGD_Echospeed ▪ Sgdperf ▪ iDrive Pro ▪ 3D FSE ▪ Smart Prep 2000 Upgrade
M3033NB 3.0T ScanTools 12.0	▪ cine ▪ epli ▪ fastgrass ▪ fse2 ▪ pcvi ▪ gtof ▪ acgdplus ▪ flairapi ▪ special ▪ ssfse ▪ threeplane ▪ e3dtof ▪ fsexl ▪ bloodsupp ▪ fastcine ▪ idrive ▪ t2bhold ▪ ssfsemrcp ▪ t1bhold ▪ ushorttr ▪ ssfsexl ▪ smartprep ▪ sgd ▪ sgdperf ▪ idrivePro ▪ FR3DFSE ▪ sprep99	▪ Cine ▪ Echo Planar Imaging ▪ Fast Gradient Echo ▪ Fast Spin Echo & FLAIR ▪ Phase Contrast Vascular Imaging ▪ Time of Flight ▪ ACGD PLUS ▪ FLAIR EPI ▪ SPECIAL ▪ SSFSE ▪ Three Plane Localizer ▪ 3D Time of Flight ▪ FSE-XL ▪ Blood Suppression ▪ Fast Cine ▪ iDrive ▪ T2 Breath Hold ▪ SSFSE MRCP/Body Pak ▪ T1 Breath Hold ▪ Ultrashort TR/TE ▪ MRCP3 ▪ Smart Prep ▪ SGD_Echospeed ▪ Sgdperf ▪ iDrive Pro ▪ 3D FSE ▪ Smart Prep 2000 Upgrade

Option Catalog	Options Activated	Customer Protocol
M3335KH 1.5T ScanTools 14.0	▪ cine ▪ epli ▪ fastgrass ▪ fse2 ▪ pcvi ▪ gtof ▪ acgdplus ▪ flairepi ▪ special ▪ ssfse ▪ threeplane ▪ e3dtof ▪ fsexl ▪ bloodsupp ▪ fastcine ▪ idrive ▪ t2bhold ▪ ssfsemrcp ▪ t1bhold ▪ ushorttr ▪ ssfsexl ▪ DynR1 ▪ smartprep ▪ asset ▪ sgd ▪ sgdperf ▪ idrivePro ▪ FR3DFSE ▪ sprep99 ▪ lava ▪ dwepi	▪ Cine ▪ Echo Planar Imaging ▪ Fast Gradient Echo ▪ Fast Spin Echo & FLAIR ▪ Phase Contrast Vascular Imaging ▪ Time of Flight ▪ ACGD PLUS ▪ FLAIR EPI ▪ SPECIAL ▪ SSFSE ▪ Three Plane Localizer ▪ 3D Time of Flight ▪ FSE-XL ▪ Blood Suppression ▪ Fast Cine ▪ iDrive ▪ T2 Breath Hold ▪ SSFSE MRCP/Body Pak ▪ T1 Breath Hold ▪ Ultrashort TR/TE ▪ MRCP3 ▪ Dynamic R1 ▪ Smart Prep ▪ Asset ▪ SGD_Echospeed ▪ Sgdperf ▪ iDrive Pro ▪ 3D FSE ▪ Smart Prep 2000 Upgrade ▪ LAVA ▪ Diffusion Weighted epi

Option Catalog	Options Activated	Customer Protocol
M3335KG 3.0T ScanTools 14.0	▪ cine ▪ epli ▪ fastgrass ▪ fse2 ▪ pcvi ▪ gtof ▪ acgdplus ▪ flairepi ▪ special ▪ ssfse ▪ threeplane ▪ e3dtof ▪ fsexl ▪ bloodsupp ▪ fastcine ▪ idrive ▪ t2bhold ▪ ssfsemrcp ▪ t1bhold ▪ ushorttr ▪ ssfsexl ▪ smartprep ▪ asset ▪ sgd ▪ sgdperf ▪ idrivePro ▪ FR3DFSE ▪ sprep99 ▪ lava ▪ dwepi	▪ Cine ▪ Echo Planar Imaging ▪ Fast Gradient Echo ▪ Fast Spin Echo & FLAIR ▪ Phase Contrast Vascular Imaging ▪ Time of Flight ▪ ACGD PLUS ▪ FLAIR EPI ▪ SPECIAL ▪ SSFSE ▪ Three Plane Localizer ▪ 3D Time of Flight ▪ FSE-XL ▪ Blood Suppression ▪ Fast Cine ▪ iDrive ▪ T2 Breath Hold ▪ SSFSE MRCP/Body Pak ▪ T1 Breath Hold ▪ Ultrashort TR/TE ▪ MRCP3 ▪ Smart Prep ▪ Asset ▪ SGD_Echospeed ▪ Sgdperf ▪ iDrive Pro ▪ 3D FSE ▪ Smart Prep 2000 Upgrade ▪ LAVA ▪ Diffusion Weighted epi
M3333TL 3T Excite ScanTools	▪ cine ▪ epli ▪ fastgrass ▪ fse2 ▪ gtof ▪ sgd ▪ flairepi ▪ special ▪ ssfse ▪ smartprep ▪ threeplane ▪ e3dtof ▪ fsexl ▪ sgdperf ▪ idrive ▪ ushorttr ▪ t2bhold ▪ t1bhold ▪ ssfsexl ▪ sprep99 ▪ FR3DFSE ▪ idrivePro ▪ acgdplus	▪ Cine ▪ Echo Planer Imaging ▪ Fast Gradient Echo ▪ Fast Spin Echo & Flair ▪ Time of Flight ▪ SGD_Echospeed ▪ Flair EPI ▪ SPECIAL ▪ SSFSE ▪ Smart Prep ▪ Three Plane Localizer ▪ 3D Time of Flight ▪ FSE-XL ▪ Sgdperf ▪ iDrive ▪ Ultra Short TR ▪ T2 Breath Hold ▪ T1 Breath Hold ▪ MRCP3 ▪ Smart Prep 2000 Upgrade ▪ 3D FSE ▪ iDrive Pro ▪ ACGD Plus
M3333TP Excite 3T iDrive Pro Plus	▪ Rtca	▪ iDrive Pro Plus

Option Catalog	Options Activated	Customer Protocol
M3333WG 3T Probe-P	▪ probep	▪ Probe-P
M3333WH 3T Probe Press and Steam	▪ spectro	▪ Probe Press and Steam
M3333WJ 3T Probe 2D-CSI	▪ probe99	▪ probe99
M3333WK 3T Probe 3D-CSI	▪ probe3d	▪ Probe3D
M3333WL 3T Research Proton Spectro	▪ Mns ▪ probep	▪ Multi-Nuclear Spectro ▪ Probe-P
M3333TJ 3T Excite 4KW MNS Collector	▪ Mns ▪ sage ▪ probep	▪ Multi-Nuclear Spect ▪ Spectro Analysis GE ▪ Probe-P
M3333TK 3T Excite 8KW MNS Collector	▪ Mns ▪ sage ▪ probep	▪ Multi-Nuclear Spect ▪ Spectro Analysis GE ▪ Probe-P
M3333TR G3 Head MNS Upgrade	▪ Mns ▪ sage ▪ probep	▪ Multi-Nuclear Spect ▪ Spectro Analysis GE ▪ Probe-P
M3333TS G3 Body MNS Upgrade	▪ Mns ▪ sage ▪ probep	▪ Multi-Nuclear Spect ▪ Spectro Analysis GE ▪ Probe-P
M3335PW SIGNA-HDx 32 Channel	▪ channel32	▪ 32-Channel
M3335PY SIGNA-HDx 1.5T HDX 16 to 32 CH UG		

Option Catalog	Options Activated	Customer Protocol
M3335RD SIGNA-HDx EX to HDX 1.5 8CH HW UG F M3335RD SIGNA-HDx EX to HDX 1.5 8CH HW UG M M3335RF SIGNA-HDx EX to HDX 1.5 8CH HW T-R M3335RG SIGNA-HDx EX to HDX 1.5 NONLCC 8CH UG M3335RJ SIGNA-HDx PRE-EX to HDX 8 CH HW UG	▪ acgdplus ▪ bloodsupp ▪ cine ▪ DynR1 ▪ e3dtof ▪ epli ▪ fastcine ▪ fastgrass ▪ flairapi ▪ FR3DFSE ▪ fse2 ▪ fsexl ▪ gtof ▪ idrive ▪ idrivePro ▪ pcvi ▪ sgd ▪ sgdperf ▪ smartprep ▪ special ▪ spre99 ▪ ssfse ▪ ssfsemrcp ▪ ssfsexl ▪ t1bhold ▪ t2bhold ▪ threeplane ▪ ushorttr ▪ lava ▪ dwepi ▪ asset ▪ asset91 ▪ merge2d ▪ bravo ▪ channel8	▪ ACGD Plus ▪ Blood Suppression ▪ Cine ▪ Dynamic R1 ▪ 3D Time of Flight ▪ Echo Planer Imaging ▪ Fast Cine ▪ Fast Gradient Echo ▪ Flair EPI ▪ 3D FSE ▪ Fast Spin Echo & Flair ▪ FSE-XL ▪ Time of Flight ▪ iDrive ▪ iDrive Pro ▪ Phase Contrast Vascular Imaging ▪ SGD_EchoSpeed ▪ SGDPerf ▪ Smart Prep ▪ SPECIAL ▪ Smart Prep 2000 upgrade ▪ SSFSE ▪ SSFSE MRCP/Body Pak ▪ MRCP3 ▪ T1 Breath Hold ▪ T2 Breath Hold ▪ Three Plane Localizer ▪ Ultra Short TR ▪ Lava ▪ Diffusion Weighted epi ▪ Asset ▪ Asset Plus ▪ 2D Merge ▪ Bravo ▪ 8-Channel

Option Catalog	Options Activated	Customer Protocol
M2033JM Openspeed4 Scantools	▪ Epli ▪ fastgrass ▪ cine ▪ fse2 ▪ gtof ▪ pcvi ▪ proactive ▪ special ▪ ssfse ▪ threeplane ▪ e3dtof ▪ fsexl ▪ idrivePro ▪ idrive ▪ ifilter ▪ FatWater	▪ Echo Planer Imaging ▪ Fast Gradient Echo ▪ Cine ▪ Fast Spin Echo & FLAIR ▪ Time of Flight ▪ Phase Contrast Vascular Imaging ▪ Proactive Service ▪ SPECIAL ▪ SSFSE ▪ Three Plane Localizer ▪ E3DTOF ▪ FSE_XL ▪ iDrive Pro ▪ iDrive ▪ Image Filter ▪ Fat Water 2001
M2090DM Openspeed4 LSDI	▪ LSDI	▪ Line Scan Diffusion Imaging
M2090SA Openspeed4 SmartPrep/step	▪ smartprep ▪ spre99	▪ Smart Prep ▪ Smart Step
M2090EE Openspeed4 DW-EPI	▪ dwepi	▪ EPI Diffusion Imaging
M20382BP Ovation4 SW Tools	▪ lfilter	▪ Image Filter
M20372BP Ovation4 Scantools	▪ fastgrass ▪ cine ▪ fse2 ▪ gtof ▪ proactive ▪ ssfse ▪ threeplane ▪ e3dtof ▪ fsexl ▪ fastcine	▪ Fast Gradient Echo ▪ Cine ▪ Fast Spin Echo & FLAIR ▪ Time of Flight ▪ Proactive Service ▪ SSFSE ▪ Three Plane Localizer ▪ E3DTOF ▪ FSE_XL ▪ Fast Cine
M20302BP Ovation4 SmartPrep/step	▪ smartprep ▪ spre99	▪ Smart Prep ▪ Smart Step
M20332BP Ovation4 EPI	▪ epli ▪ dwepi	▪ Echo Planer Imaging ▪ EPI Diffusion Imaging
M20342BP Ovation4 Fiesta 2D/3D	▪ fiesta2d ▪ fiesta3d ▪ fiestac	▪ 2D Fiesta ▪ 3D Fiesta ▪ Phased Cycled Fiesta
M20352BP Ovation4 Fat/Water Separation	▪ FatWater	▪ Fat Water 2001
M20362BP Ovation4 Rapid Frame	▪ RapidFrame	▪ Rapid Frame
M20322BP Ovation4 LSDI	▪ LSDI	▪ Line Scan Diffusion Imaging
M20392BP Ovation4 Idrive Pro Plus	▪ idrivePro ▪ idrvie ▪ rtca	▪ iDrive Pro ▪ iDrive ▪ iDrive Pro Plus
M3335PW SIGNA-HDx 32 Channel		
M3335PY SIGNA-HDx 1.5T HDX 16 to 32 Ch UG	▪ channel32	▪ 32-Channel

5- GLOSSARY

Bloodsupp - Blood Suppression

Cine - Images usually of moving anatomy scanned and arrange for playback rapidly in movie format.

ConnectPro (HIS/RIS) - Provides DICOM modality work list connectivity to hospital or radiology information system.

DW EPI - Diffusion Weighted Echo Planer Imaging

eOM – electronic Order Management, the Oracle order management system.

Fast cine - Fast gated multi-phase images of the heart that show the heart as it contracts and expands during cardiac cycle in a movie format.

FTP - Short for *File Transfer Protocol*, the protocol used to exchange files between computer systems.

IVI - Interactive Vascular Imaging.

iDRIVE Pro - Interactive imaging uses a fast GRE pulse sequence and a new user interface to acquire continuous 2D images at a high frame rate and provide on-image navigation. Not available on systems with a single i860 AP.

Multi-Nuclear Spectroscopy - Provides access to additional spectroscopy PSD that may be used in proton mode.

MS/MA - Multi-Slice/Multi-Angle Oblique

Sgdperf - SGD enhancements, which increase, slew rates and duty cycles for higher resolution thin slice imaging.

SmartPrep - Automated triggering for contrast enhanced MRA

SmartPrep2000 - Provides SmartPrep automated table motion

Spectroscopy/PROBE - Single voxel proton spectroscopy with PROBE-S and PROBE-P

Probe2000 Upgrade - Multi voxel proton spectroscopy (PROBE - SI)

SPECIAL - Spectral / Spatial Fat Suppression. Chemical saturation using an inversion pulse (180° pulse). This key will always suppress just fat and cannot be changed to suppress water.

SSFSE - Single Shot Fast Spin Echo

Tagging - The placement of "tags" in an MR image is used to assess myocardial wall

functions visually.