

Condition	Reference	Effectivity
Removal of some magnet enclosure components may be required to complete Magnet Electrical Checks. Refer to either of the following documents and the appropriate Service Methods documentation (available through the Common Documentation Library at www.gehealthcare.com or your local GE Healthcare Service Representative) before continuing with this procedure if cover removal is required.	Chapter 5, Upper and Side Wide Open Enclosure Removal	-
-	-	-

5.4 Procedure

1. Engage the Shim Lead Assembly in conformance with Shim Lead Engagement and Disengagement.
2. Locate the connector pins using [Table 4-5](#), [Table 4-6](#), [Table 4-7](#) and [Illustration 4-20](#).

NOTE: Wiring is shown in [Chapter 6, Magnet System Wiring](#).

3. Use a digital meter to measure the resistance across the identified connector pins.
4. Measure the resistances and record in the Measured column: Main Coil resistances in [Table 4-5](#), Super Conducting (S/C) Shim Coil resistances measured on the Shim Lead in [Table 4-6](#) and S/C Shim Coil resistances measured on the Shim Lead Pigtail Cable in [Table 4-7](#).
5. Report all Measured readings outside the acceptable range (typical values listed in Tables 1, 2 and 3) to the Regional Magnet & Cryogenics (MAC) Team Leader for further investigation before ramping the magnet.

Illustration 4-20: Connector Pinouts

Shim Lead (Main Extension Lead) Housing

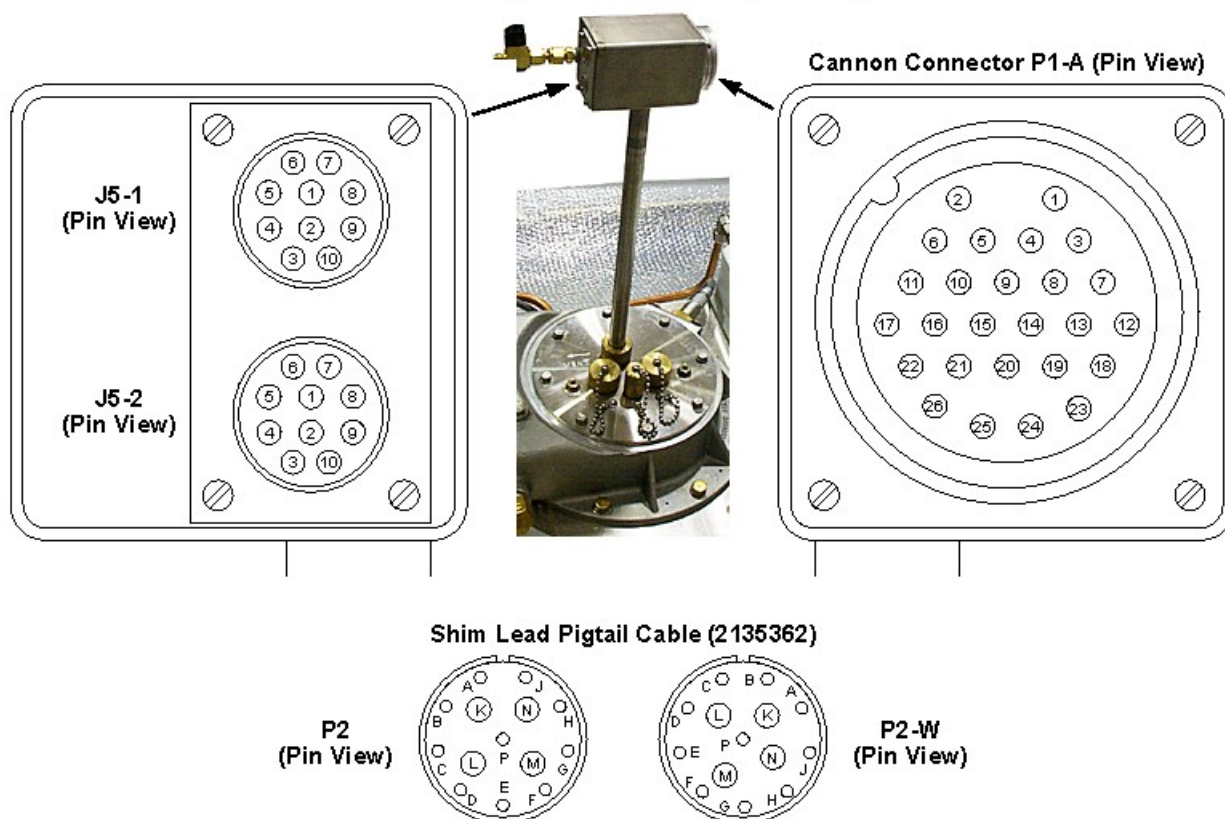


Table 4-5: Main Coil Resistance Checks, Cold (4.2K)

Connector	Pins	Resistance (Ohms)	
	+, -	Typical	Measured
Main Coil Power Lugs	+, -	< 6	
J5-1 on Shim Lead	9, 10	< 6	

Table 4-6: S/C Shim Coil Resistance Checks, Cold (4.2K) - Measured on Shim Lead

Function		Connector	Pins	Resistance (Ohms)		
Shim PS Set-ting	Shim Coil Cir-cuit		+, -	Typical		Measured
				R, RA	RB	
AX1	Z1	P1-A	1, 19	0.3 - 0.5	0.04 - 0.80	
AX2	Z2		2, 20			
AX3	Z3		3, 21			
AX4	Z4		4, 22			
AX5	Z5		5, 23			
AX6	Z6		6, 24			
T1-1	C31		13, 23			
T1-2	C11+		16, 19			
T1-3	C22+		14, 21			
T1-4	C11-		17, 20			
T1-5	C22-		15, 22			
T1-6	C33		18, 24			
T2-1	S31		11, 23			
T2-2	S11+		9, 19			
T2-3	S22+		7, 21			
T2-4	S11-		10, 20			
T2-5	S22-		8, 22			
T2-6	S33		12, 24			
S/C Switch Heaters Main Switch		J5-1 or J5-2	1, 2	22 - 27		
Axial Shims			5, 6	27 - 33 or 22 - 28*	18 - 22	
Transverse 1			7, 6	8 - 10	12 - 18	
Transverse 2			8, 6			

* New switch configuration incorporated around R216 and Q063.

* New switch configuration incorporated around R216 and Q063.

Table 4-7: S/C Shim Coil Resistance Checks, Cold (4.2K) - Measured on Shim Lead Pigtail Cable

Function		Connector	Pins	Resistance (Ohms)		
Shim PS Set- ting	Shim Coil Cir- cuit		+, -	Typical		Measured
				R, RA	RB	
AX1	Z1	P2	A*, (K,L)	0.3 - 0.6	0.04 - 0.80	
AX2	Z2	P2-W				
AX3	Z3	P2	B, (M,N)			
AX4	Z4	P2-W				
AX5	Z5	P2	C, P			
AX6	Z6	P2-W				
T1-1	C31	P2	D, P			
T1-2	C11+		F, (K,L)			
T1-3	C22+		G, (M,N)			
T1-4	C11-	P2-W	F, (K,L)			
T1-5	C22-		G, (M,N)			
T1-6	C33		D, P			
T2-1	S31	P2	E, P			
T2-2	S11+		H, (K,L)			
T2-3	S22+		J, (M,N)			
T2-4	S11-	P2-W	H, (K,L)			
T2-5	S22-		J, (M,N)			
T2-6	S33		E, P			
* New switch configuration incorporated around R216 and Q063.						

5.5 Finalization

No finalization steps.