

[new...](#) | [new revision](#) | [page status](#) | [copy](#) | [move](#) | [delete](#)

welcome!
ct tech pubs web page
ct water cooler
ct tech pubs updates
surveys
9800/hla
highspeed family (us)
ct/i
lightspeed family
yms family
all ct products
index
old ct emails
communiques pending
approval
customize
members

search:

advanced search
news: daily | weekly

go offline | chat | notify |
print | tutorial | help

home> all ct products

Table Encoder Cable Interference

christopher hardiman (christopher c hardiman), 05/31/2002 - 08:15 pm



PROPRIETARY WARNING

This information is Property of GE - For GE Service Personnel Only.
Disclosure to third parties prohibited.

Communique Number

MQUE-5ANGVW

Equipment

All HiSpeed and LightSpeed Products Using the Zeus Table

Problem

The table longitudinal encoder cable can be snagged near the gantry end of the table by customer-introduced obstructions. This often results in damage to the encoder assembly through "sudden spool unwind" and the cable being knocked off of the spool.

Solution

A new cover has been designed to protect the cable and add some additional clearance for the spool. This was to be introduced on the LightSpeed Plus systems; however, 214 tables were manufactured using the old cover. It is recommended that the old cover (46-297897P1) be replaced with this new cover (2261675). It was the result of a Six Sigma project and was found to dramatically reduce the frequency of these types of failures. The cover is compatible with and can be installed on all flavors of the Zeus table. The cost is minor and it is recommended that you order the new cover to replace your old cover on your next pm. Drawings of the two covers are included below for reference.

You can download any of the files below by dragging them to your Windows desktop.
To open a file for reading, double-click it.




Author

Christopher C Hardiman

Revision

05/31/2002 08:15 PM


[new...](#) | [new revision](#) | [page status](#) | [copy](#) | [move](#) | [delete](#)

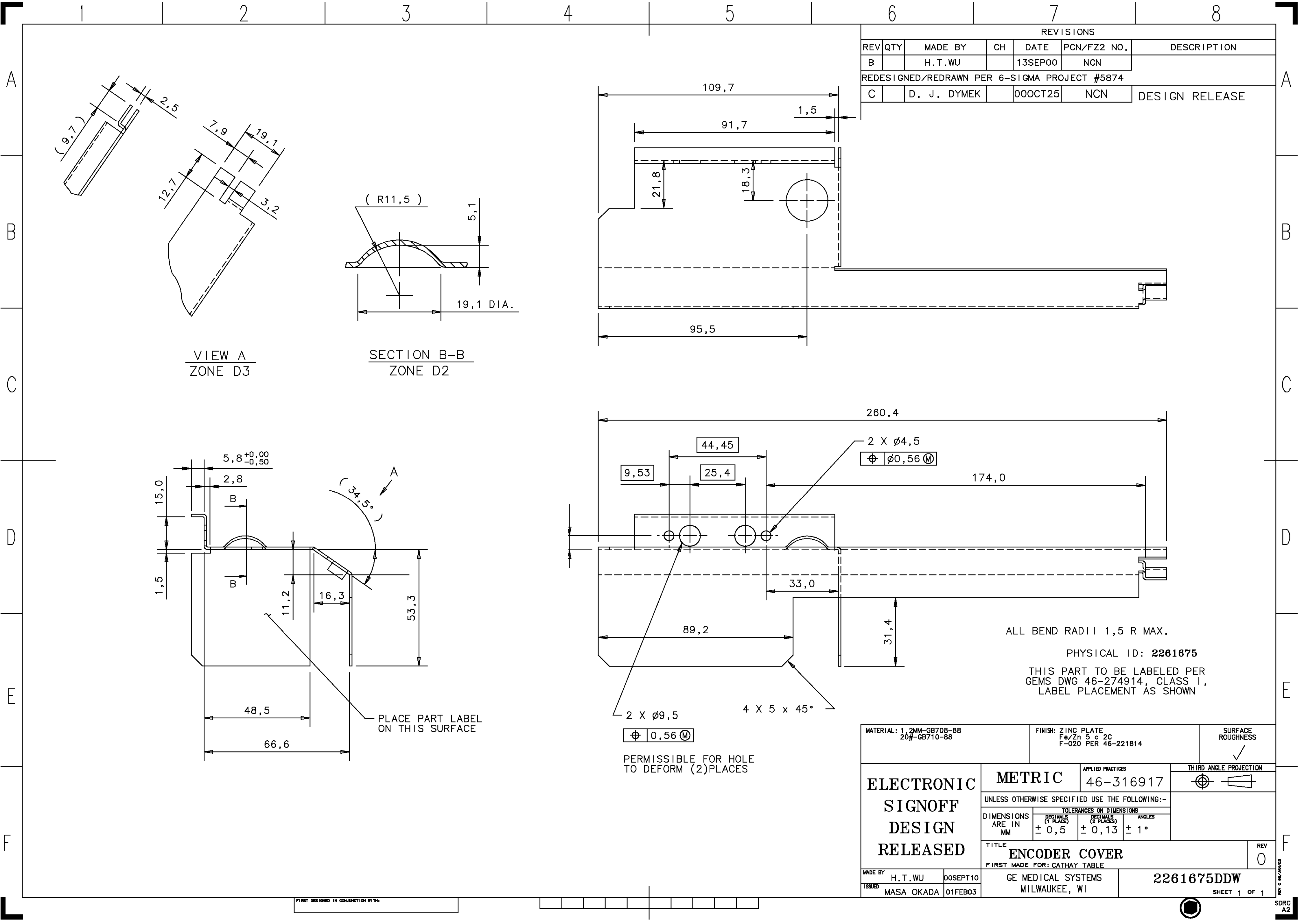

46-297897
DRAWING NO. 46-297897
CONT. ON SHEET 1

GE Medical Systems

46-297897
CONT. ON SHEET 1

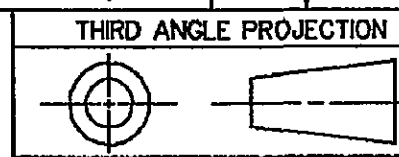

THIRD ANGLE PROJECTION

UNLESS OTHERWISE SPECIFIED USE THE FOLLOWING:-

APPLIED PRACTICES

46-170700

SURFACES

✓

TOLERANCES ON DIMENSIONS

DECIMALS
2 PLACES
± .02




DECIMALS
3 PLACES
±

ANGLES
± 5°

TITLE

ENCODER COVER

FIRST MADE FOR ZEUS TABLE

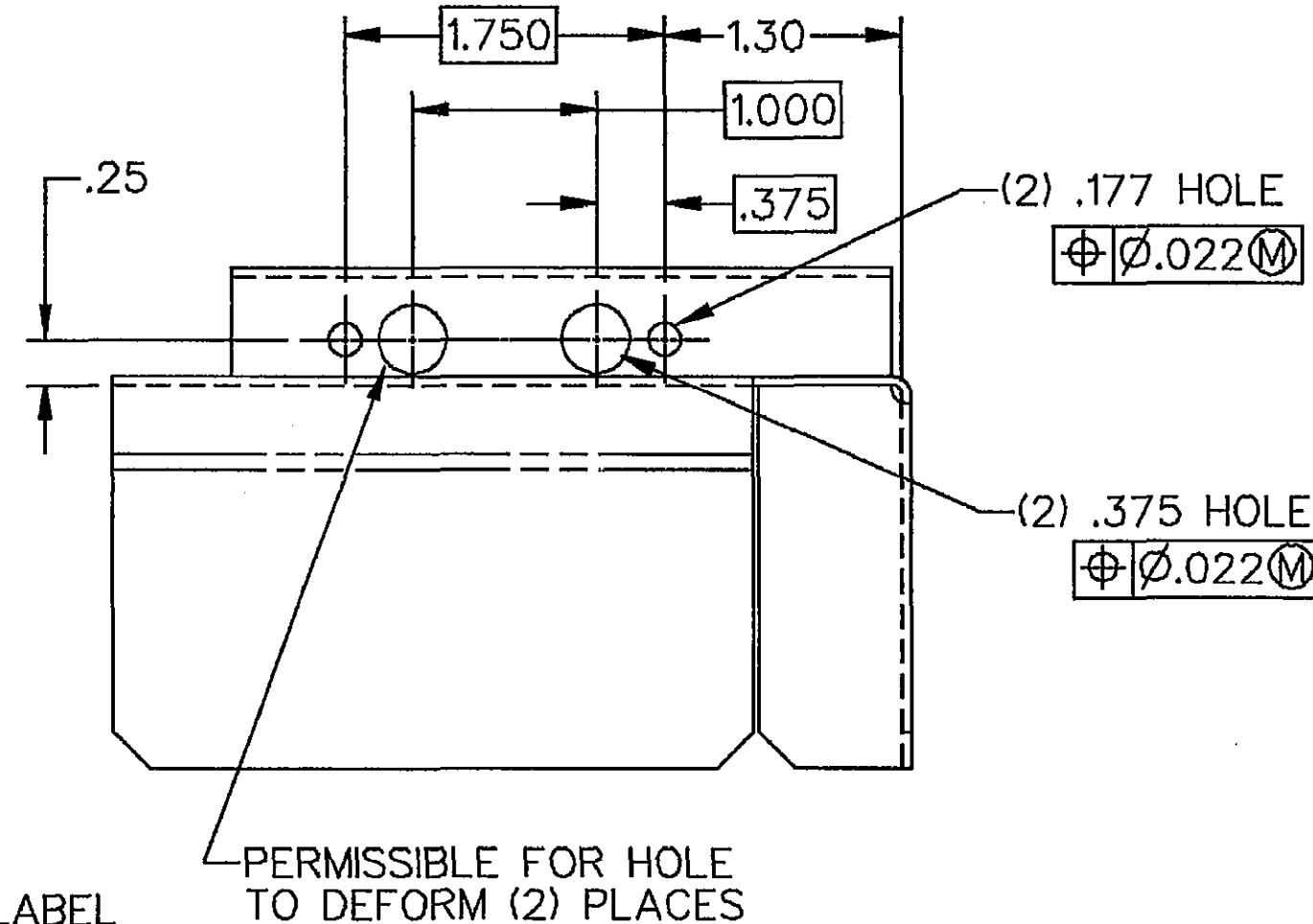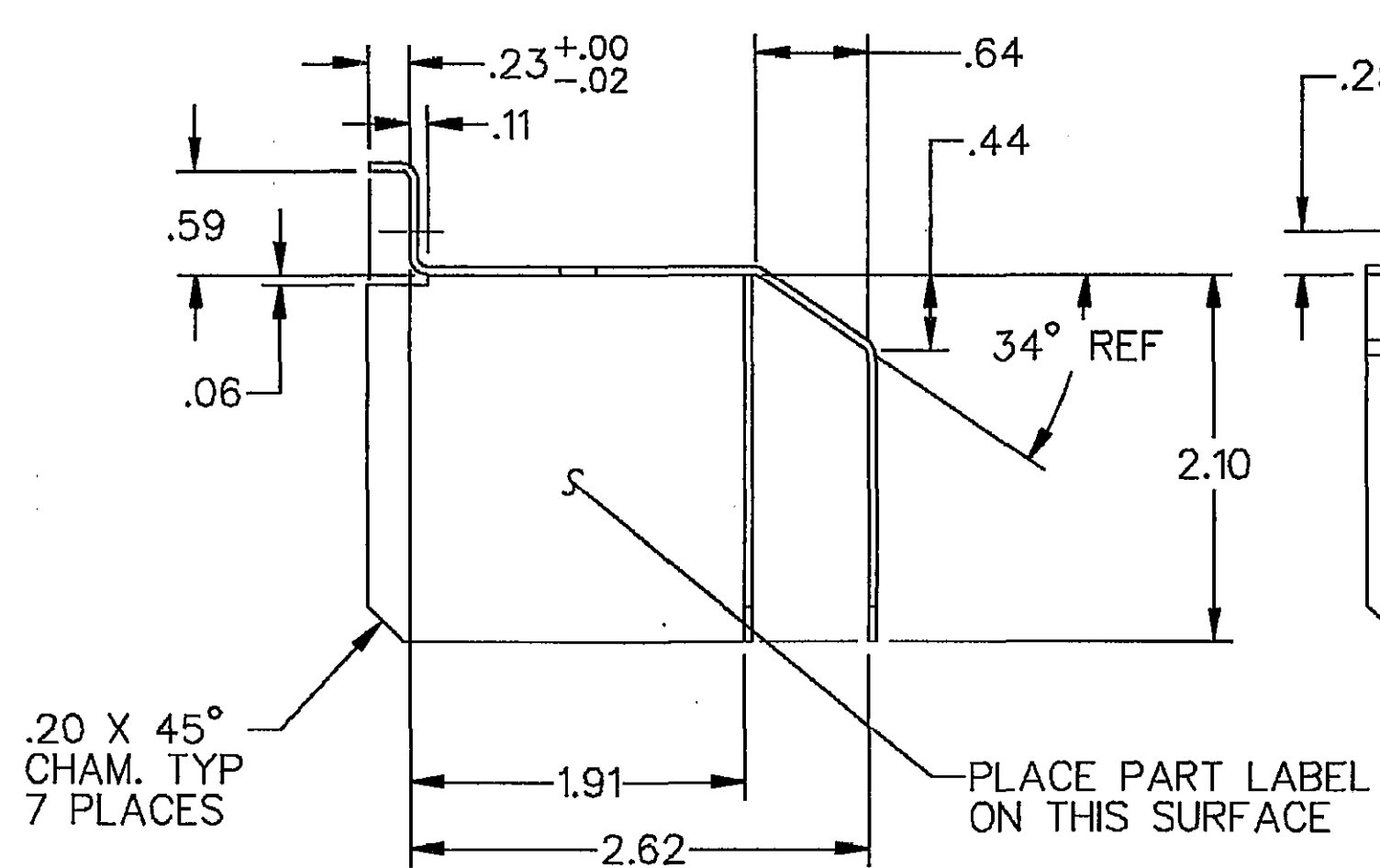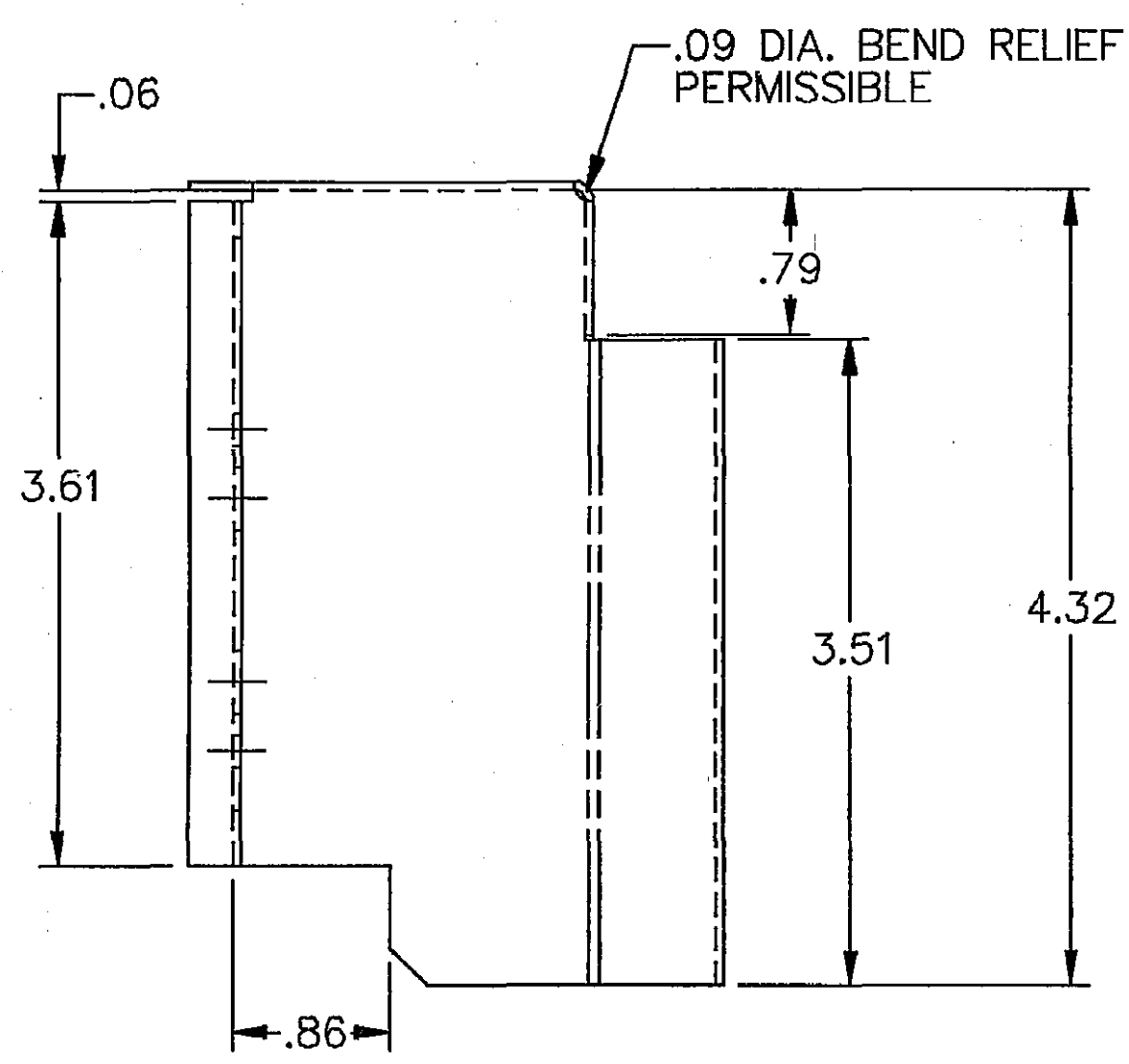

ALL BEND RADII .06 R MAX.

- 1 MATL: C.R. SHEET STEEL, AISI 1010
ASTM A366, A619, AND A620
(REF. GE EMPIS B8A4)
.048 THICK (18 GAGE)
FINISH: BRIGHT ZINC ELECTROPLATE
(.0002-.0004 THICK)
(REF. GE EMPIS F70B5A)
F-020 PER GEMS DWG 46-221814

THIS PART TO BE LABELED PER
GEMS DWG 46-274914, CLASS I,
LABEL PLACEMENT AS SHOWN

APPROVALS		
	NAME	DATE
TECH	DMK	6/23/92
DES	JFR	6/22/92
ENG	C VARA	6/22/92
MFG	GEG	6/22/92

MADE BY
JOE ROBINETTE 92APR28
6-22-92

REVISIONS			PRINTS TO	
B	JOE ROBINETTE 92JUN19	NCN	500	
REMOVED WELD, .06, .16 & .11, REVISED PART LABEL, .79-WAS-.73, ADDED .09 DIA. BEND RELIEF PERM. .25 & HOLE DEFORM PERM. REMOVED .34, .23-WAS-.24				

GE MEDICAL SYSTEMS
MILWAUKEE, WISCONSIN

DIV OR
DEPT.
LOCATION

46-297897
CONT. ON SHEET 1