

[new...](#) | [edit](#) | [page status](#) | [copy](#) | [move](#) | [delete](#)

welcome!
ct tech pubs web page
ct water cooler
ct tech pubs updates
surveys
9800/hla
highspeed family (us)
ct/i
lightspeed family
yms family
all ct products
[index](#)
old ct emails
communiques pending
approval
communique change
requests
customize
members

search:

advanced search
news: daily | weekly

go offline | chat | notify |
print | tutorial | help

home> index

Brush Block Guide Pins

christopher hardiman (christopher c hardiman), 09/20/2002 - 10:28 am

**Priority**

3 - Troubleshooting Information

PROPRIETARY WARNING

This information is Property of GE - For GE Service Personnel Only.
Disclosure to third parties prohibited.

Communique Number

MQUE-5E6JEV

Equipment

LightSpeed Plus, Ultra & 16

Problem

The brushes on the slipring brush block are susceptible to breakage if lateral force is applied during assembly back to the gantry.

Solution

A guide pin, P/N: 2358034, has been specified and tested to aid in the placement of the brush blocks. This pin is a FRU and two should be ordered for each LightSpeed Plus, Ultra and 16 system an FE services. These pins will be added to the Installation Support Kit shipped with each system in the near future. For more information, please refer to the below document.

You can download any of the files below by dragging them to your Windows desktop.
To open a file for reading, double-click it.

**Author**

Christopher C Hardiman

Revision

09/20/2002 10:28 AM

[new...](#) | [edit](#) | [page status](#) | [copy](#) | [move](#) | [delete](#)



GE Medical Systems

Technical Publication

**Direction 2352788-100
Revision 0**

LightSpeed Plus, Ultra & 16, Brush Block Guide Pins

Copyright 2002 by General Electric Company



Revision History

REV	DATE	REASON FOR CHANGE
Rev 0	08/08/2002	Initial Release

List of effective pages

Page	Rev.	Page	Rev.	Page	Rev.
Title Page	0	1	0	Rear	-
Rev/LOEP	0				

Problem:

The brushes on the slipring brush block are susceptible to breakage if lateral force is applied during assembly back to the gantry.

Solution:

Included in the Installation Support kits for LightSpeed Plus, Ultra and 16 systems are two M6 x 47.5 mm studs (see Figure 1, p/n: 2358034) to be used as guide pins. Use of these “pins” will aid in the placement of the brush block and prevent damage to the brush tips.



Figure 1. M6 x 47.5 mm Stud

Remove the top and bottom mounting bolts and insert these threaded rods in their place (see Figure 2). Remove the other two bolts and slide the brush block out using the studs as guides. Perform the necessary cleaning, or removal of the block, and reverse the procedure to re-install the brush block. The guide pins should be removed and the M6 mounting bolts installed in their place upon reinstallation.

Do not leave the M6 studs on the system.

Be sure to torque the (4) 6 mm cap screws to 7.9 N-m (70 in-lbs or 6 ft-lbs).

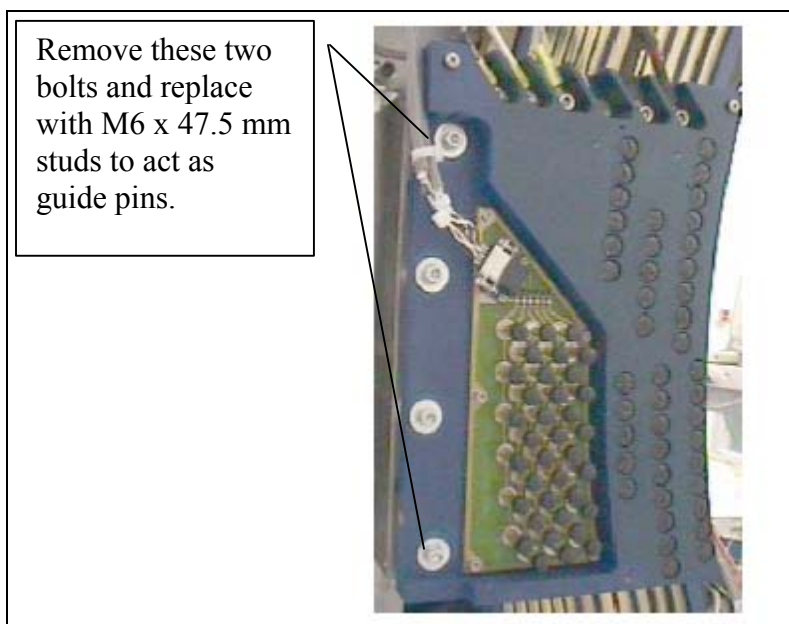


Figure 2. Brush Block

Keep these pins on site for future use during cleaning!



GE Medical Systems

*CT SERVICE ENGINEERING, W-1200
P.O. BOX 414
MILWAUKEE WISCONSIN 53201-0414 U.S.A.
Fax: (262) 312-7434*