

GE HealthCare

Technical Bulletin



Technical Bulletin Unique Identifier	DOC2848327
Title	Service Note for the impact of removal of ICC facility bypass valve in 1.5T SIGNA Prime & 1.5T SIGNA Victor/ Star Air/ Aviator Air
Subject	Service Note for the impact of removal of ICC facility bypass valve in 1.5T SIGNA Prime & 1.5T SIGNA Victor/ Star Air/ Aviator Air
Applicability (e.g., affected product, process)	1.5T SIGNA Victor 1.5T SIGNA Aviator Air 1.5T SIGNA Star Air 1.5T SIGNA Prime
Date	07-Aug-2023
Relevant Product Business (select one)	☒ Imaging ☐ Patient Care Solutions ☐ Ultrasound ☐ Other: _____
Originating Modality	MR
Details	To avoid the problem of cooling water backflow, Option 1: Switch facility water to cooled cryogenic compressor as Emergency Backup Water Supply method is now removed.
Tracking Number	ECR2356897
Affected Service Publications	1.5T SIGNA Victor/Star Air/Aviator Air Service Manual Class A: 5957356-8EN Rev.1; Class C: 5957355-8EN Rev.1; Class M: 5957347-8EN Rev.1 1.5T SIGNA Prime Service Manual Class A: 5807565-2EN Rev.4; Class C: 5807569-2EN Rev.4; Class M: 5807568-2EN Rev.3
Resolution	Keep Option2: Switch city water to cooled cryogen compressor as the only Emergency Backup Water Supply backup method for ICC cabinet (5842254). Latest ICC cabinet removed the 2 water valves BV11 and BV12, and the rest sides have been blocked.
Technical Bulletin Owner	Qin Yu (SSO:503343456)

(This Service Note is classified as Class A)

*** Printed copies are uncontrolled unless otherwise identified ***

Before using this document, consult MyWorkshop for the latest revision.
GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.

GE HealthCare

Technical Bulletin



Effectivity:

1.5T SIGNA Prime

1.5T SIGNA Victor / 1.5T SIGNA Aviator Air / 1.5T SIGNA Star Air

Background:

There are 2 options for F-50SH Emergency Backup Water Supply:

Option 1: Switch facility water to cooled cryogen compressor

Option 2: Switch city water to cooled cryogen compressor


While applying **Option 1: Switch facility water to cooled cryogenic compressor** for Emergency Backup Water Supply it may cause the cooling water pour back.

To avoid this problem, **remove Option 1: Switch facility water to cooled cryogenic compressor** as Emergency Backup Water Supply method, keep **Option2: Switch city water to cooled cryogen compressor** as the only Emergency Backup Water Supply backup method for ICC cabinet (5842254).

The design of the latest ICC cabinet (5842254 Rev.3) removed the 2 water valves BV11 and BV12, and the rest sides of them have been blocked.

This Service Note is for declaring the impact of removal of ICC facility bypass valve in 1.5T SIGNA Prime & 1.5T SIGNA Victor/ Star Air/ Aviator Air.

Figure 1. Removed Valve BV11 & BV12

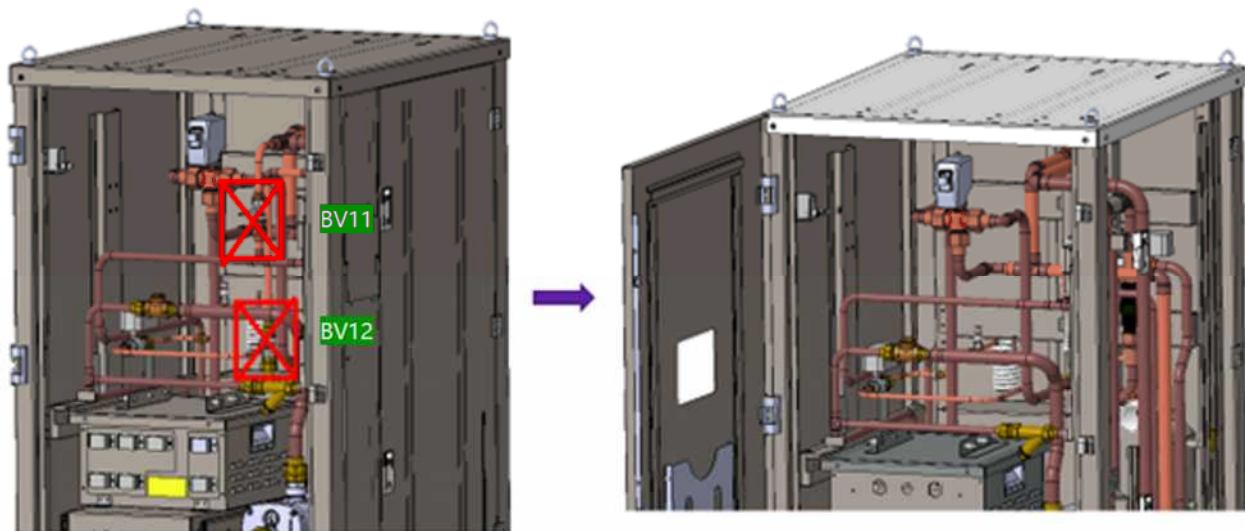

GE HealthCare

Technical Bulletin



Impact:

Service Manual:

Functional Checks/ System / Power on sequencing / Integrated Cooling Cabinet (ICC) power on procedure

Before power on:

For the ICC cabinet (5842254 Rev.3), it is no longer necessary to check whether BV11 & BV12 valves are closed.

Service Manual:

Troubleshooting / Cooling / F-50SH Emergency Backup Water Supply

There are 2 options for F-50SH Emergency Backup Water Supply:


Option 1: Switch facility water to cooled cryogen compressor

Option 2: Switch city water to cooled cryogen compressor

Removed Option 1: Switch facility water to cooled cryogen compressor.

Option 2: Switch city water to cooled cryogen compressor is the **only way as F50-SH Emergency Backup.**

Figure 2. Removed Valve BV11 & BV12

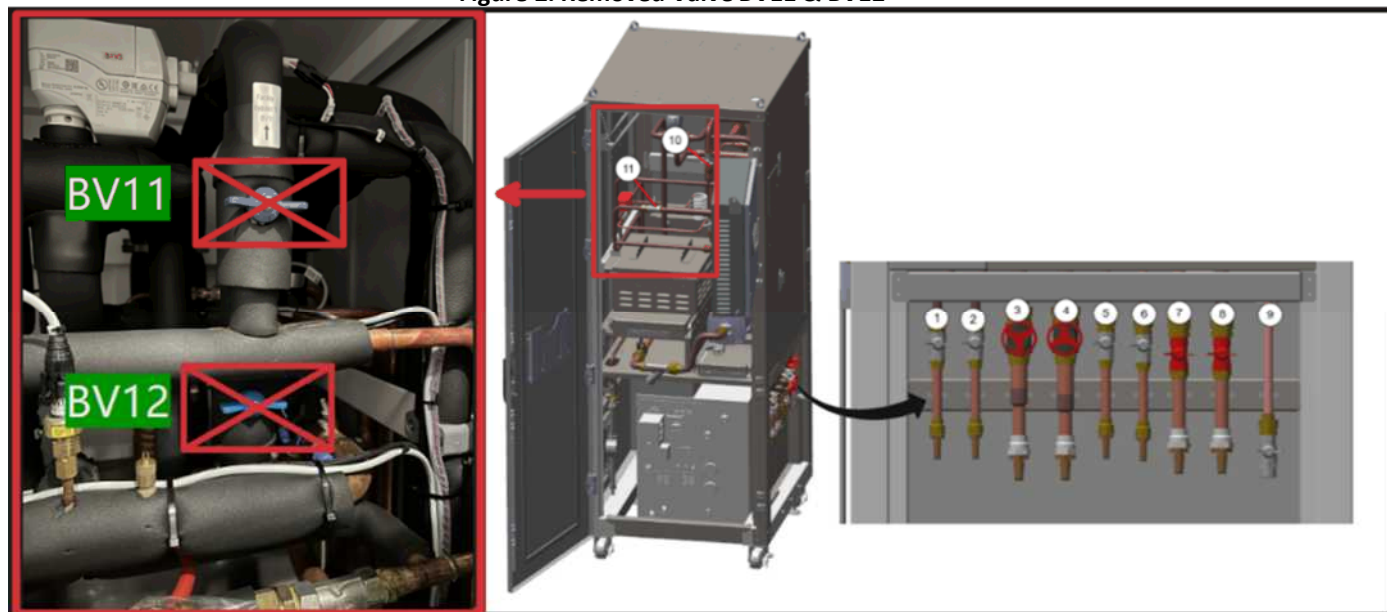

*** Printed copies are uncontrolled unless otherwise identified ***

Before using this document, consult MyWorkshop for the latest revision.
GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.

GE HealthCare

Technical Bulletin



Table 1. ICC Valve

Item	Description
1	From City Water BV09
2	To City Drain BV10
3	From Facility BV03 (Out- side Facility Water Con- trol Valve)
4	To Facility BV04 (Outside Facility Control Valve)
5	To ISC BV32
6	From ISC BV31
7	From GC BV28
8	To GC BV27
9	Tank Drain
10	To Facility BV11 (Inside Facility Water Control Valve) Add: For ICC cabinet (5842254 Rev.3) there is no BV11 valve
11	From Facility BV12 (Inside Facility Water Control Valve) Add: For ICC cabinet (5842254 Rev.3) there is no BV12 valve

The latest ICC cabinet (5842254 Rev.3) also removed the Emergency value label description on the lower right-side panel of the ICC cabinet and behind the ICC control box.

Figure 3. Removed Label in the latest ICC cabinet

*** Printed copies are uncontrolled unless otherwise identified ***
 Before using this document, consult MyWorkshop for the latest revision.
 GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.

GE HealthCare Technical Bulletin


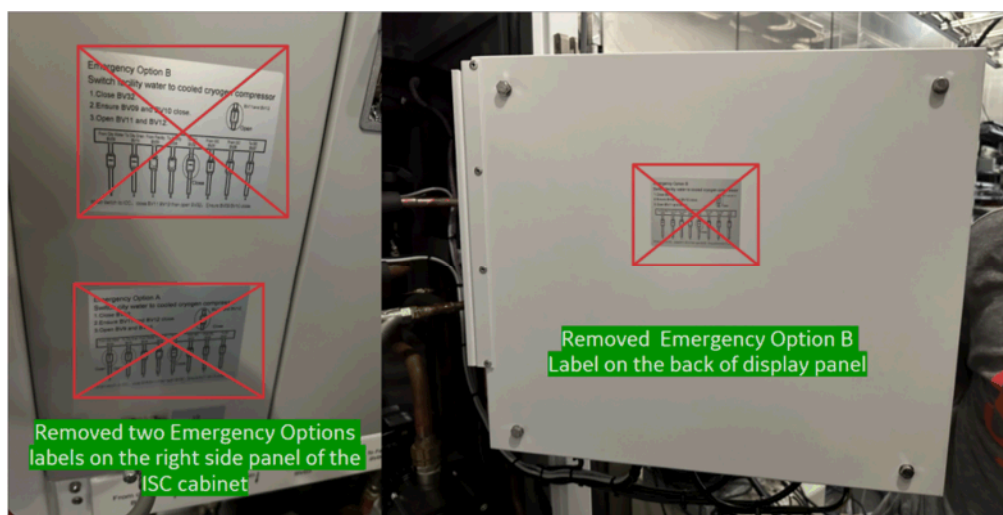


Figure 4. Add new operation guide label under value on the lower right-side panel

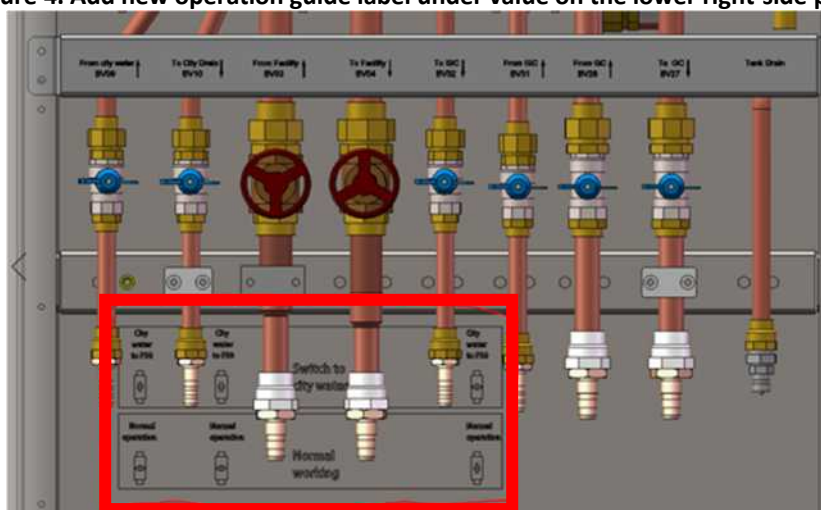

*** Printed copies are uncontrolled unless otherwise identified ***
Before using this document, consult MyWorkshop for the latest revision.
GE HealthCare Confidential


General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.
Reproduction and/or distribution is prohibited.

Released

GE HealthCare

Technical Bulletin

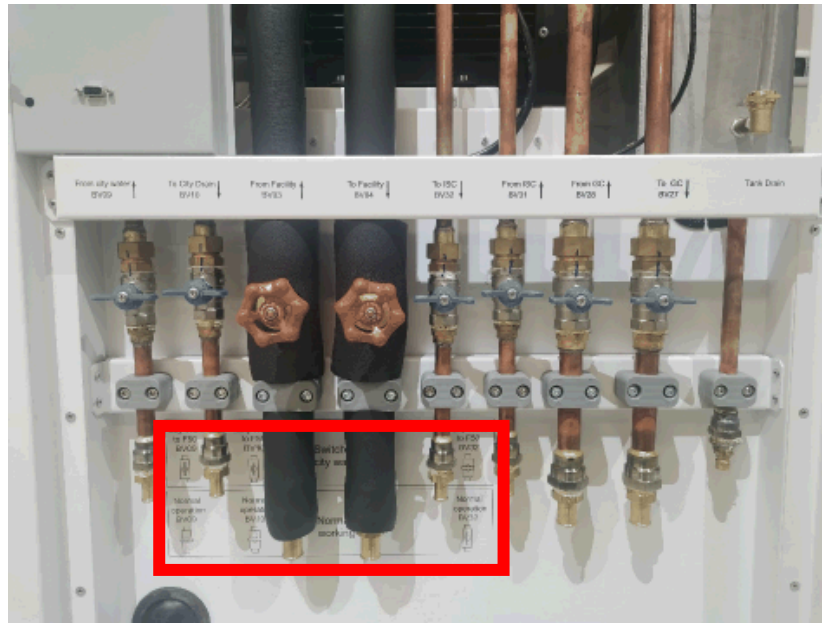

Revision	Version	Date	Author	Description
01	01	07-Aug-2023	Qin Yu (503345456)	Initial Release

End of Document

*** Printed copies are uncontrolled unless otherwise identified ***
 Before using this document, consult MyWorkshop for the latest revision.
 GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.