

GE HealthCare

Technical Bulletin



Technical Bulletin Unique Identifier	DOC2884628
Title	Service Note for content update in Section 7.2.11 of MR450w non-GEM to SIGNA Artist Upgrade Manual
Subject	Content update in Section 7.2.11 of MR450w non-GEM to SIGNA Artist Upgrade Manual
Applicability (e.g., affected product, process)	MR Optima MR450W 1.5T MR Signa Artist 1.5T
Date	29-Aug-2023
Relevant Product Business (select one)	☒ Imaging ☐ Patient Care Solutions ☐ Ultrasound ☐ Other: _____
Originating Modality	MR
Details	Service Note for content update in Section 7.2.11 of MR450w non-GEM to SIGNA Artist Upgrade Manual
Tracking Number	N/A
Affected Service Publications	MR450w non-GEM to SIGNA Artist Upgrade Manual Part #: 5791841 Rev.8 Class A
Resolution	See attached
Technical Bulletin Owner	Qin Yu (SSO: 503343456)

(This SN is classified as Class A)

*** Printed copies are uncontrolled unless otherwise identified ***

Before using this document, consult MyWorkshop for the latest revision.

GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.

GE HealthCare

Technical Bulletin



Effectivity:

MR Optima MR450W 1.5T

MR Signa Artist 1.5T


Details:

Since it was found from the site that four **(4) M10 Wedge Spacer Washers (5417957)** described in Step 1 of **Chapter 7.2.11 Install Middle Support Bracket (Dv Magnet)** does not exist in the BOM and was not found in the site, to avoid misunderstandings.

This SN updates step 1 of the MR450w to Signa Artist upgrade manual. Replace the content of the original **red text** of Step1 with the content of the following **green text**.

7.2.11 Install Middle Support Bracket (Dv Magnet)

Figure 1 (Figure 123 in manual). Middle Supports

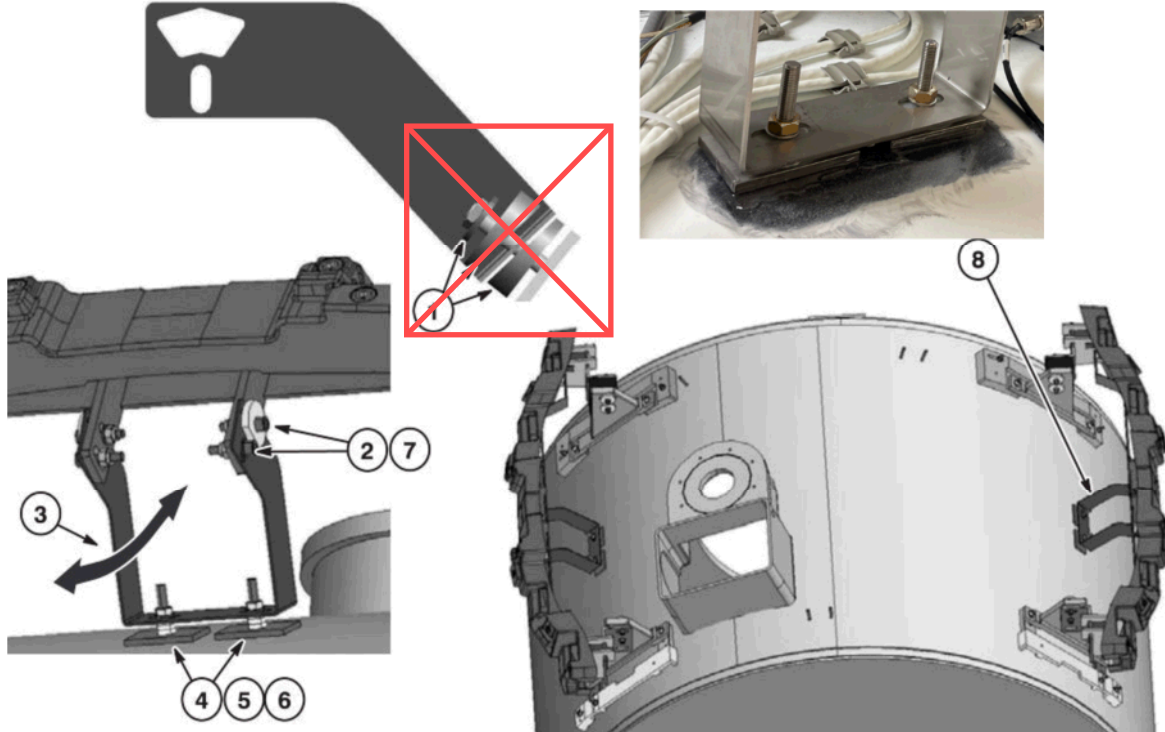

1. Assemble Middle Support Bracket (5414810), and Replacement Studs with two (2) M10 Brass Hex Nuts (3000-M10C-40), two M10 SS Flat Washers (2000-M10-10), and four (4) M10 Wedge Spacer Washers (5417957).


NOTE: Make sure the M10 Wedge Spacer Washers are aligned as shown above.

GE HealthCare

Technical Bulletin



Figure 2. The actual situation of installation after update

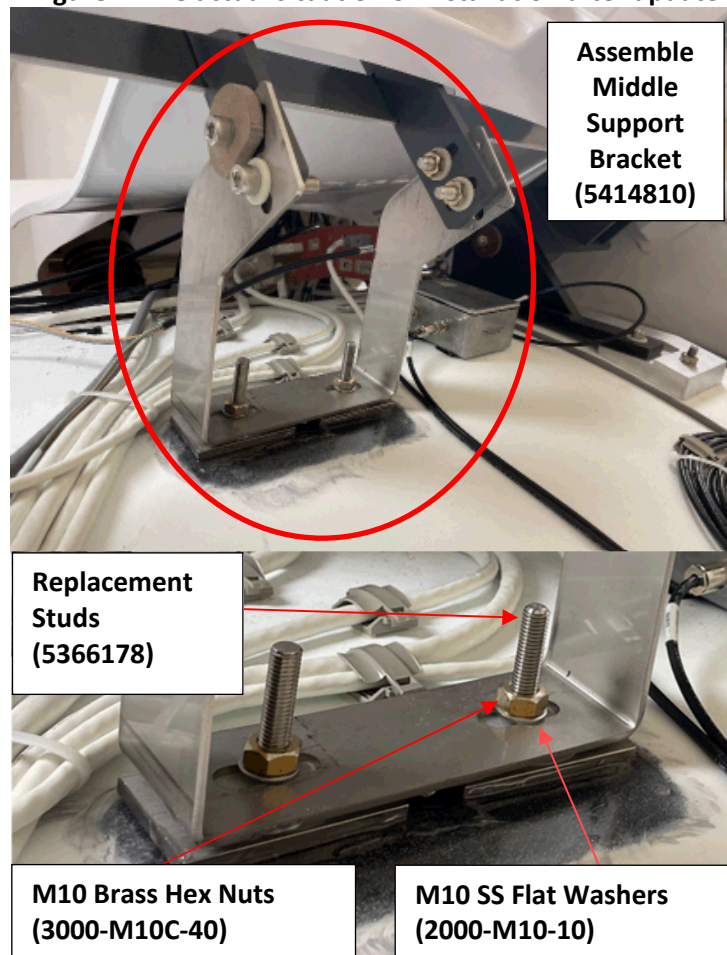

1. Assemble Middle Support Bracket (5414810), and Replacement Studs (5366178) with two (2) M10 Brass Hex Nuts (3000-M10C-40), two (2) M10 SS Flat Washers (2000-M10-10).

2. Attach Middle Support Bracket to Top Frame using two (2) M10 x 30mm HSHC Screws (1000-M10C030-33 at the top location, (2) M10 x 25mm HSHC Screws (1000-M10C025-33) at the bottom location, two (2) M10 Large Washers (5372304), six (6) M10 SS Flat Washers (2000-M10-10), and four (4) M10 Brass Hex Nuts (3000-M10C-40). Do not tighten.

3. Rotate the Middle Support Bracket to establish proper location of replacement studs on surface of magnet.

4. Mark location of replacement studs base on magnet surface.

5. Move the Middle Support Bracket to make room for removal of paint on marked surface using provided sandpaper and tack cloth. Make sure surface is clean of all excess paint particles.

*** Printed copies are uncontrolled unless otherwise identified ***

Before using this document, consult MyWorkshop for the latest revision.

GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.

GE HealthCare

Technical Bulletin



6. Align the studs back to the marked positions.

NOTE: The replacement studs will be bonded to the magnet surface **AFTER** the enclosure covers have been installed and aligned.

7. Tighten fasteners between Middle Support Bracket to Top Frame.

8. Repeat above steps for installing right side Middle Support.

Revision	Version	Date	Author	Description
01	01	29-Aug-2023	Qin Yu	Initial Release

End of Document

*** Printed copies are uncontrolled unless otherwise identified ***

Before using this document, consult MyWorkshop for the latest revision.

GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.