GE Healthcare

Signa EXCITE 4–Channel to 8–Channel Upgrade

Technical Publication

5177001
Rev 1

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation “**damage in shipment**” written on **all** copies of the freight or express bill **before** delivery is accepted or “signed for” by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

File a report with

- Call 1–800–548–3366 and use option 8.
- Contact your local service coordinator for more information on this process.

Rev. 08/15/2003

© 2006, General Electric Company
GE Medical Systems, a General Electric Company, going to market as GE Healthcare
3200 N. Grandview Boulevard
Waukesha, Wisconsin 53188
USA

www.gehealthcare.com

Table of Contents

DESCRIPTION	PAGE
Table Of Contents	1
Language Policy For Service Documentation	3
Revision History	5
SECTION 1 – Overview	7
SECTION 2 – Upgrade Procedure	9

This page left blank intentionally.

WARNING

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

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Revision History

Date	Rev	Change Description
2–14–2006	1	Initial Release

This page left blank intentionally.

SECTION 1 – OVERVIEW

1–1 INTRODUCTION

This manual provides instructions to upgrade a Signa EXCITE HD or HDx system from a 4–Channel system to 8–Channels.

This Direction covers the following catalogs:

M3000LL – 4 Channel to 8 Channel Upgrade Kit

This upgrade adds:

- A second cable track containing Runs 1210 and 1211
- The B–Bezel and harness to the LPCA

1–2 SITE READY CHECK

Pre–installation work must be completed before equipment is delivered to avoid delays and confusion.

1–3 DISCARDED MATERIAL RETURN POLICY

The equipment or material being removed by the upgrade have a potential value for service of the installed base. The disposition of such must be controlled by GE. Disposition of the removed material is per GE Medical Systems Policy and Procedures.

1–4 SOFTWARE LOAD FOR UPGRADE

Install the 8–Channel Option key.

1–5 PRODUCT LOCATOR

The Global Install Base Database (also known as Product Locator System) tracks shipment, trans–shipment, and field location of the serialized models. There are now two methods for submitting Product Locator information.

1. At this time, for **U.S. ONLY**, the preferred method of submitting information is the *FE Site Verification Web Site* at:
<http://gein2.med.ge.com/gib>

The FE Site Verification consists of three components that are available on the web from the main menu. They are:

- Install/deinstall product locator model and serial numbers
- Add/modify ship to address information
- Update CARES FE data for primary/secondary FEs

To obtain a copy of the FE training tutorial for using this website, a downloadable copy is available at the following:
Product Locator Support Central Page – <http://supportcentral.ge.com/15563>

2. One “Shipping Card” is filled out and submitted when shipped (extra cards are supplied for trans–shipments between storage and distribution points), and the “Installation Card” and extra shipment cards are attached.

Verify that serial and model number on each rating plate matches installation card numbers before removing installation card. Note that there may be one or more shipment cards and bar code labels with the installation card. These shipment cards are used to trace the transfer of serialized units between various inventory storage and distribution points until the product reaches its final installation destination. Process just the installation card and discard any extra shipment cards and labels.

1–6 FINAL COORDINATION

1. Replace all covers that have not been previously replaced.
2. Record and enter applicable data into applicable site configuration files and records.
3. Process product locator installation cards for all serialized components added to the system to include product locator card supplied with the upgrade serial number rating plate. Return to:

Product Locator File
P.O. Box 414, W–523
Milwaukee, WI 53201–0414

Note


Failure to fill out and return Product Locator Cards may result in failure of your site to receive future FMI's.



SECTION 2 – UPGRADE PROCEDURES


2-1 REMOVE LPCA COVER

① Move LPCA to rear of Bridge.


② Unfasten the four screws that attach the cover to the base plate.

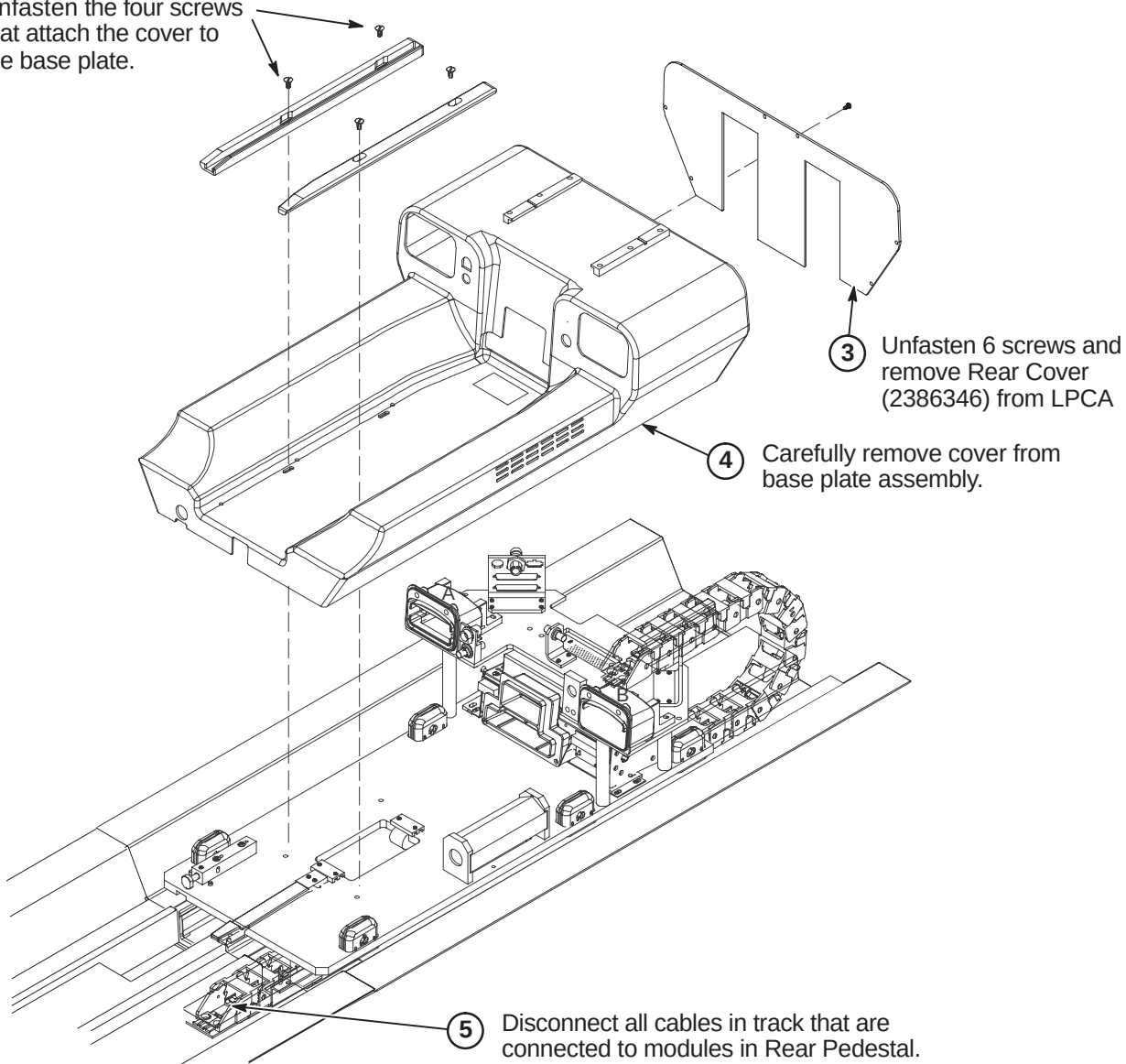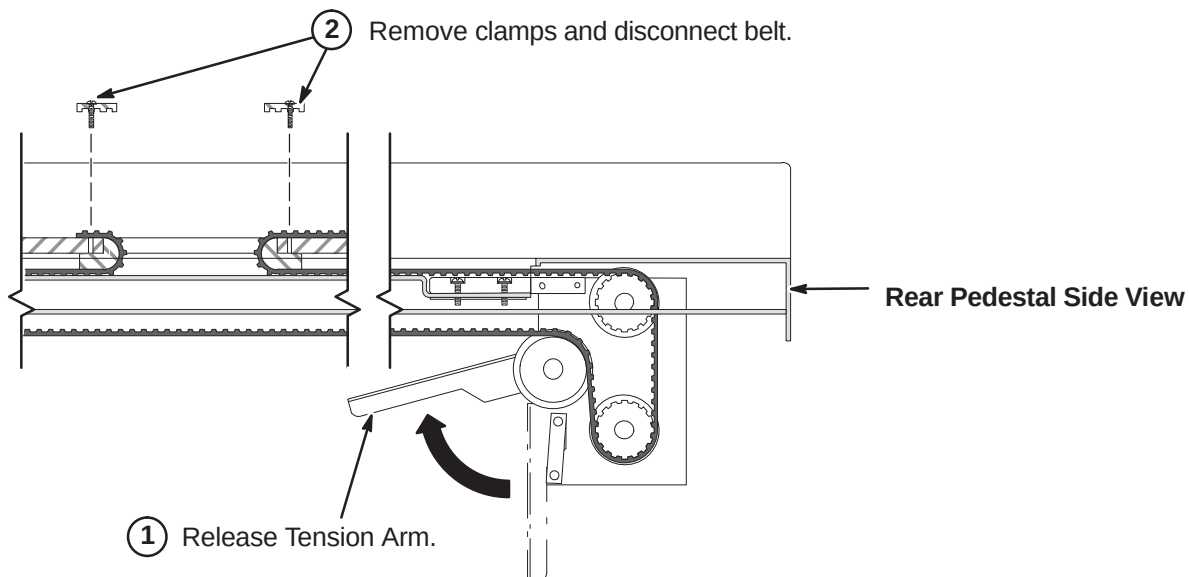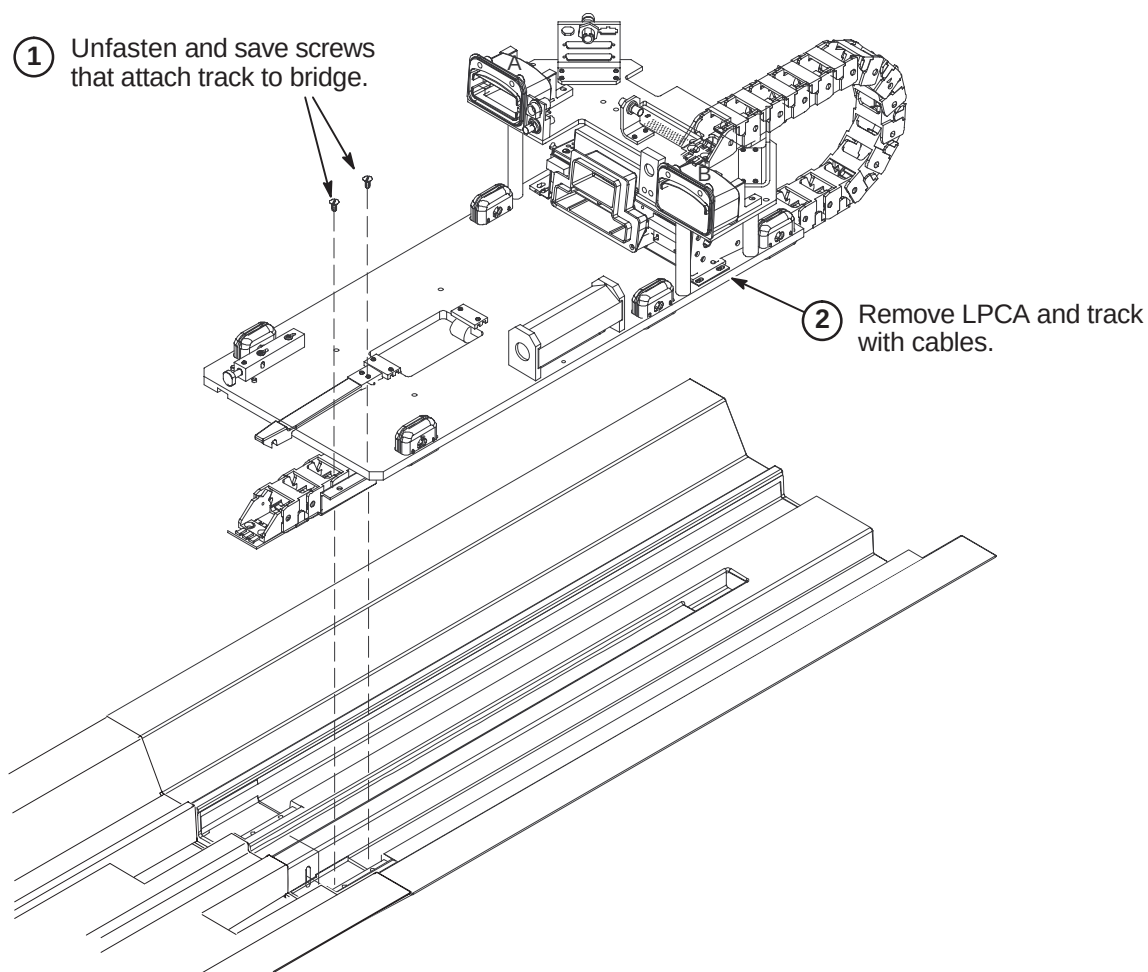
2-2 RELEASE TENSION ARM AND DISCONNECT BELT**2-2-1 Remove LPCA and Cable Track**

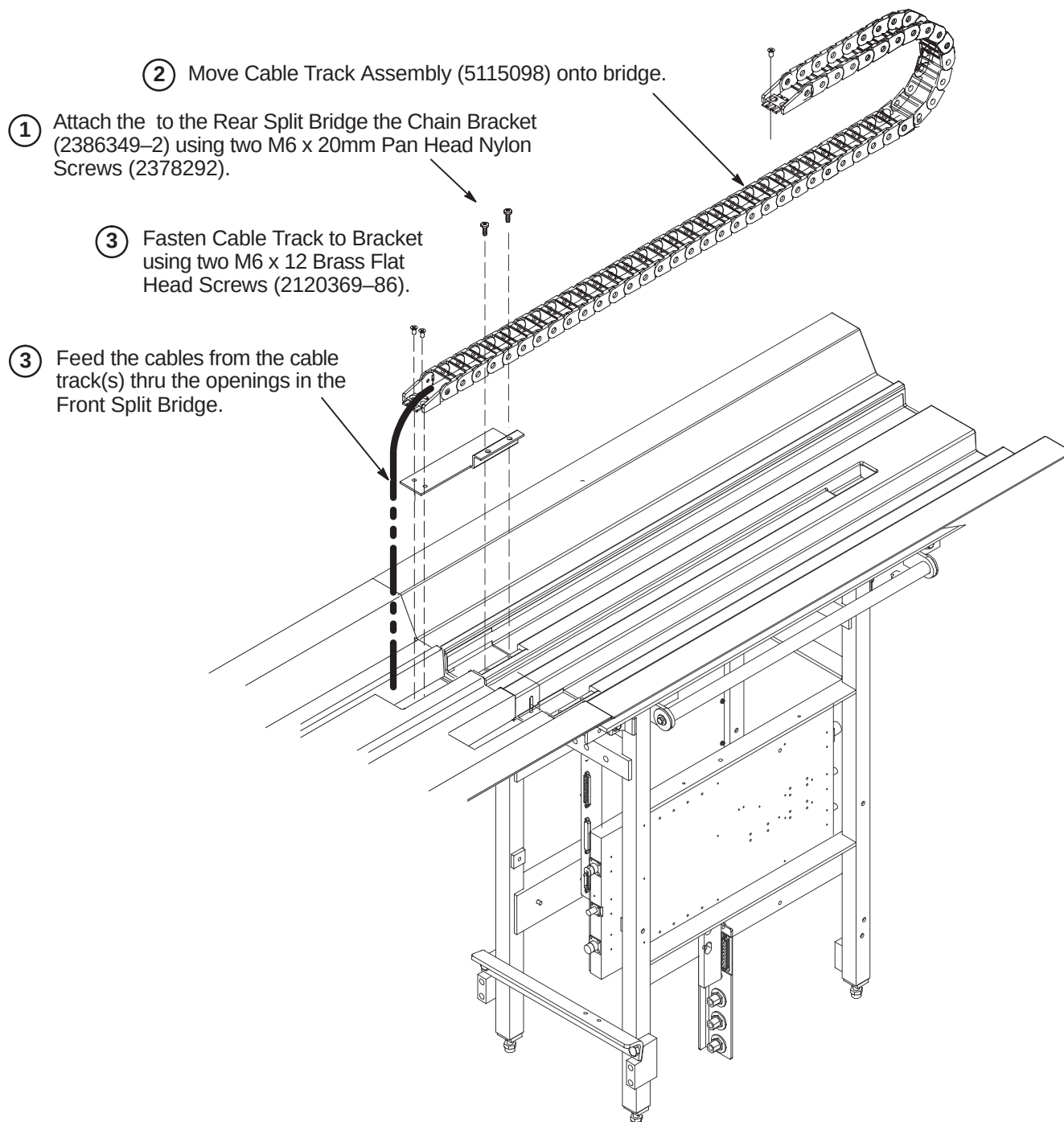
2-3 ATTACH CABLE TRACKS TO BRIDGE

2-4 ROUTE AND CONNECT RUNS 1210 AND 1211

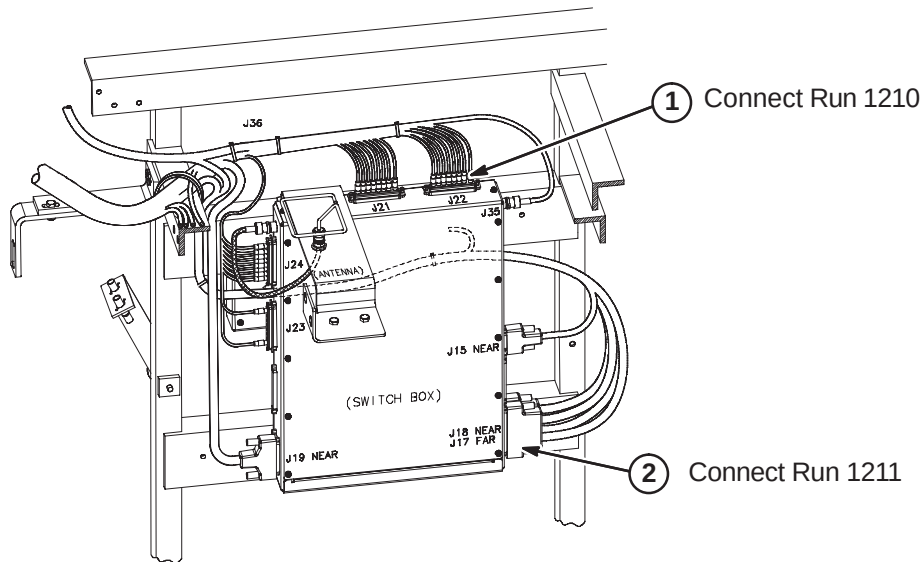
2-4-1 Install LPCA and Cable Track

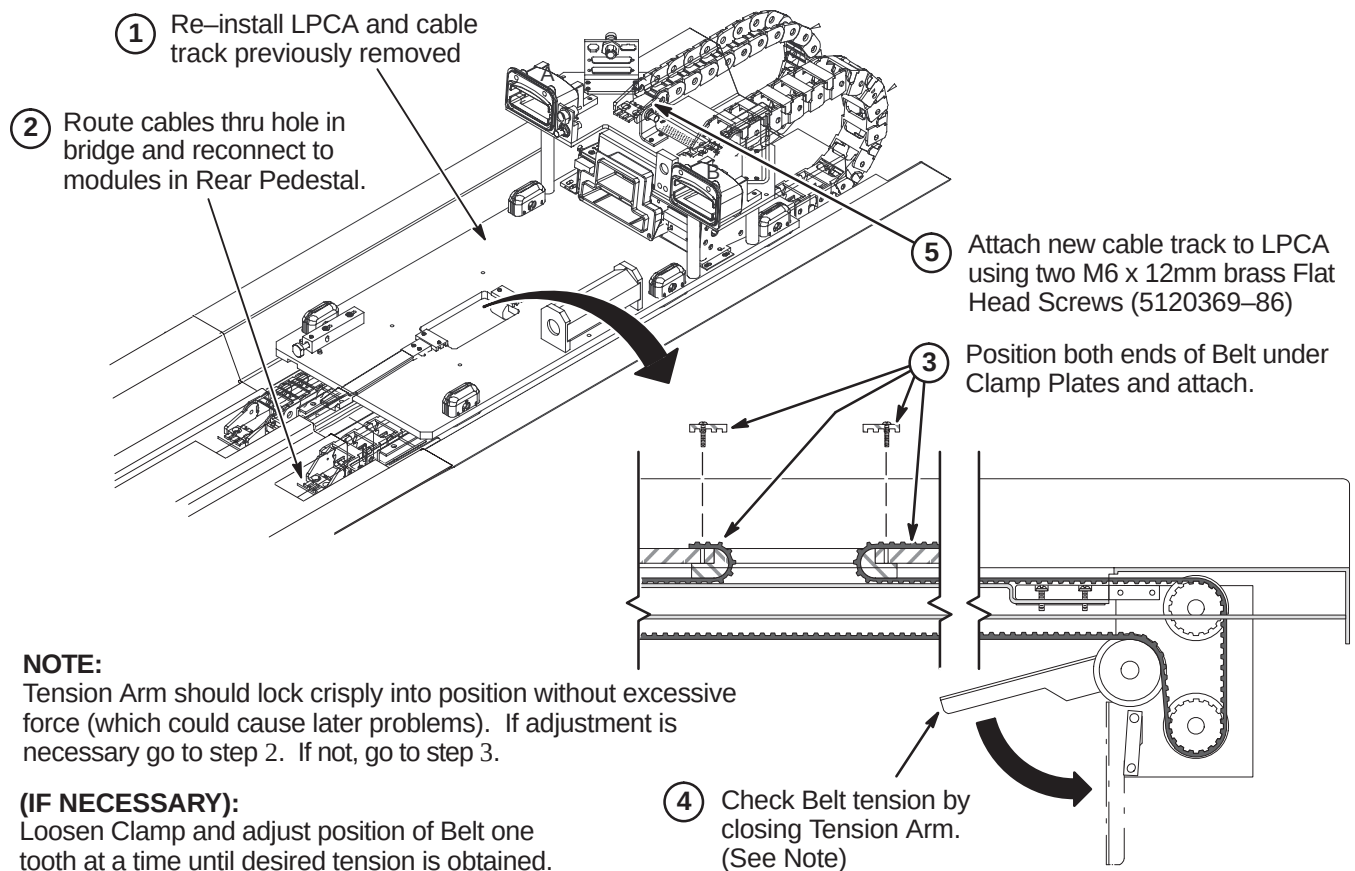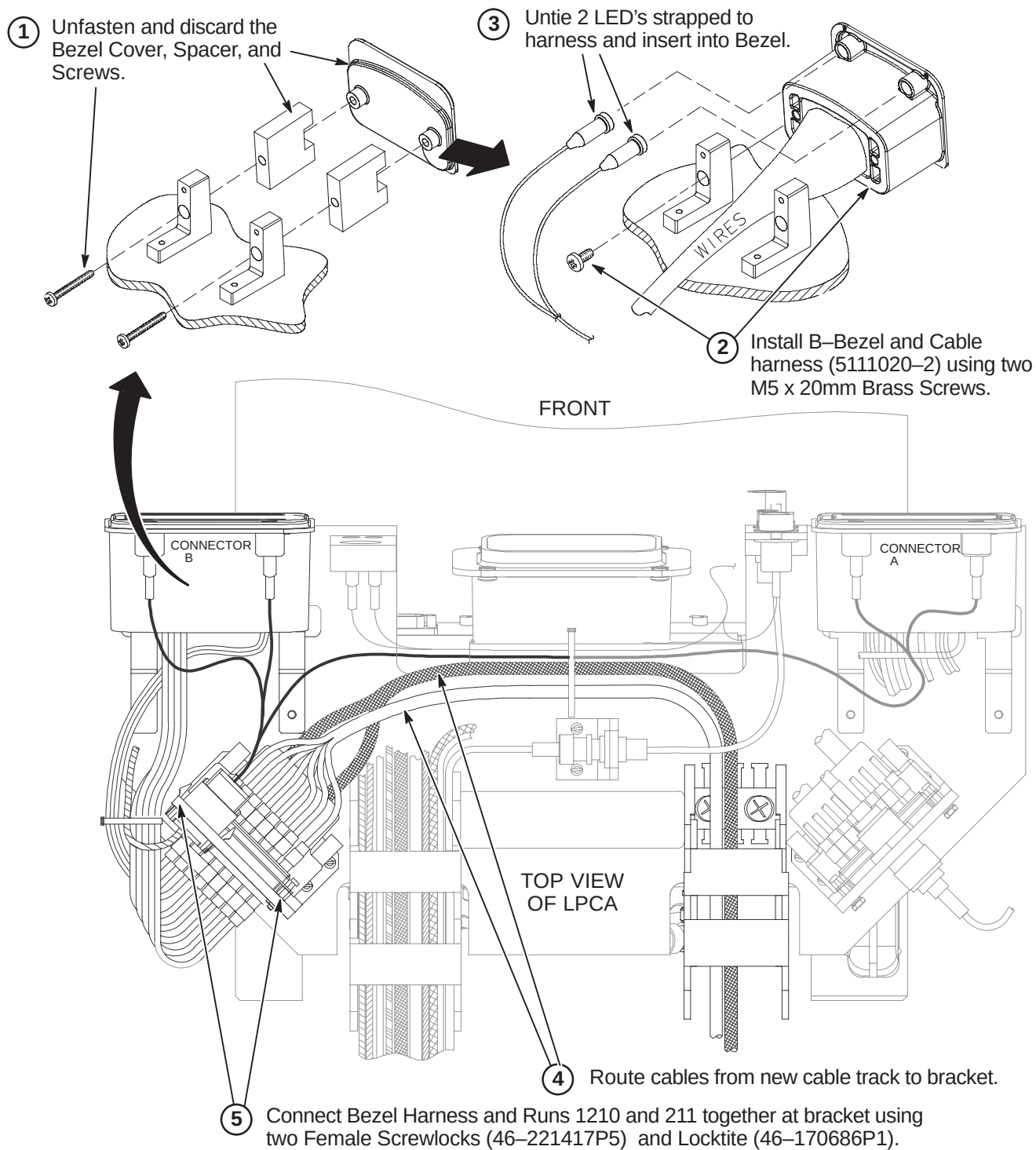2-5 INSTALL AND CONNECT B-BEZEL TO RUNS 1210 AND 1211 IN LPCA

2-6 Installation of LPCA Cover

- ③ Attach the cover to the base plate using four screws previously removed.

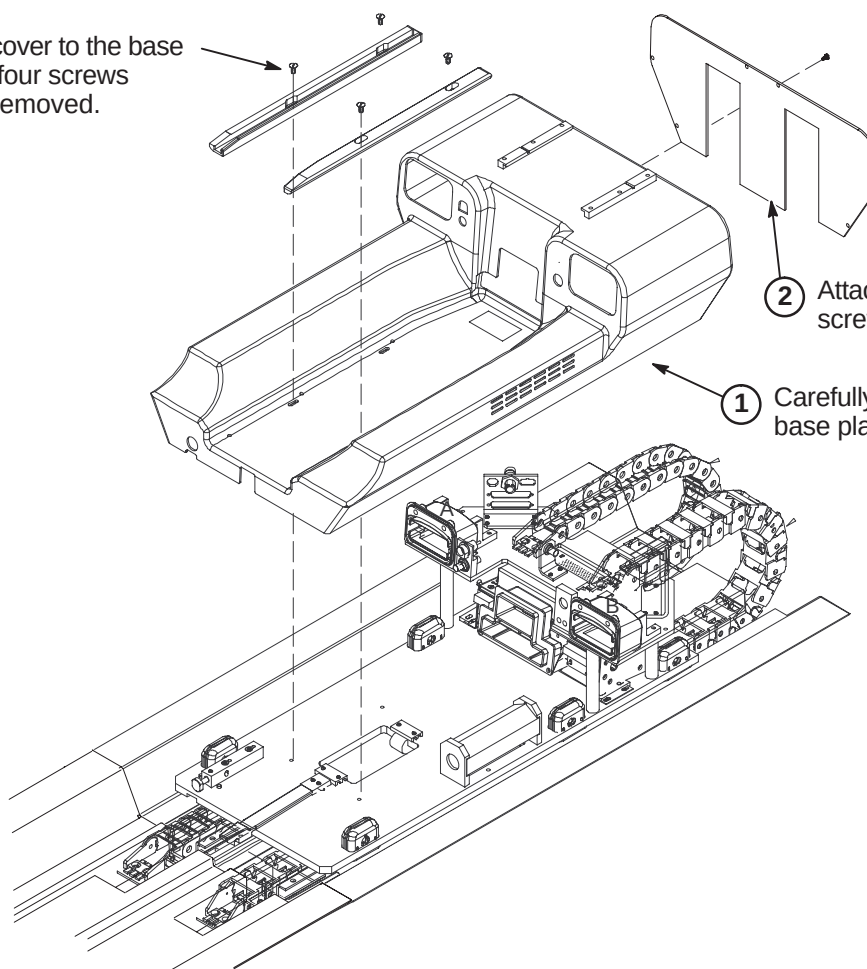- ② Attach Rear Cover using 6 screws previously removed.

- ① Carefully place cover on base plate assembly.