

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1- BODY T/R- HYBRID/COMBINER REPLACEMENT	2
1-1 Tools and Instruments Required.....	2
1-1 Removing the Body Hybrid Splitter.....	2
1-1 Installing the Body Hybrid Splitter	3
2- FINALIZATION.....	4
REVISION HISTORY.....	5

1- BODY T/R- HYBRID/COMBINER REPLACEMENT


This procedure outlines the process for replacing the 3.0T Excite Body T/R Module.

1-1 Tools and Instruments Required


1. Medium Phillips Screw Driver
2. Basic Hand tools

1-1 Removing the Body Hybrid Splitter

- 1 Make sure to End Exam.
- 2 Disable the output of the RF amplifier by using the Switch on the Exciter or Disconnecting J14 at the rear of the RF Amplifier.
- 3 Remove the two lower covers from the rear bridge assembly.
- 4 Disconnect the T/R BNC.
- 5 Disconnect the Heliax Cables.
- 6 Slide the Hybrid splitter out of its saddle.
- 7 Two People need to carry the unit from the magnet room. Make sure to stay as far away from the Magnet as possible.
- 8 Once outside of the Magnet Room, Remove the two Preamps and Notch Filter. They will be reinstalled on the new Body Hybrid/Combiner assembly.

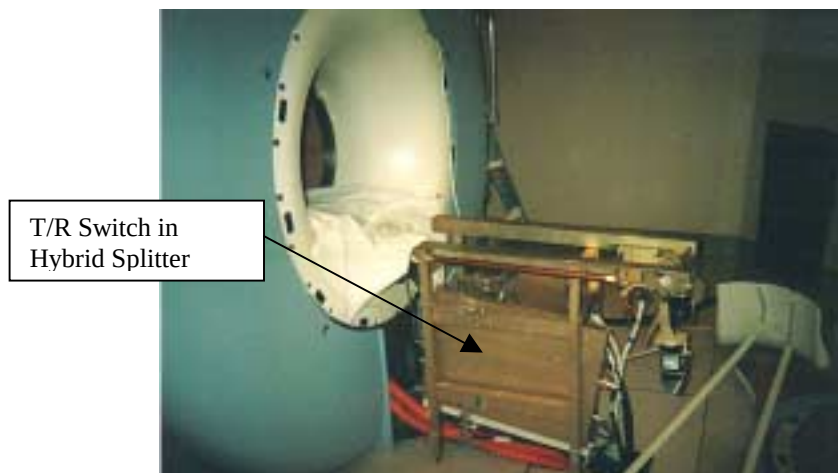


LOCATION OF HYBRID SPLITTER
ILLUSTRATION 2-1

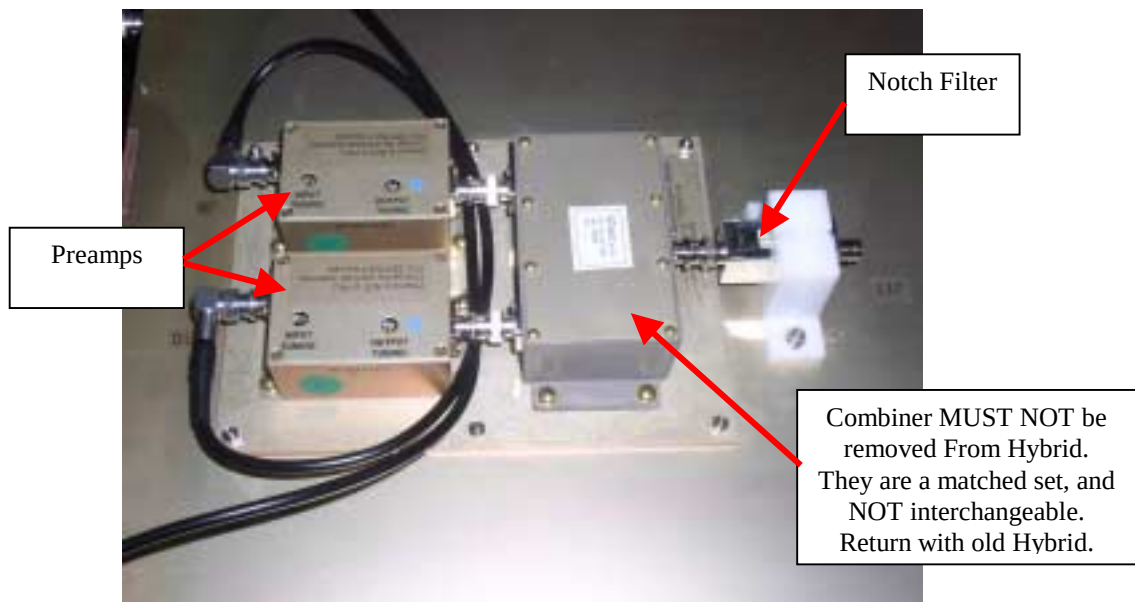

REMOVE TWO PREAMPS AND NOTCH FILTER FROM OLD UNIT
ILLUSTRATION 2-2

1-1 Installing the Body Hybrid Splitter

- 1 Attach the existing Preamps and Notch Filter to the new Hybrid Splitter/Combiner Assembly.
- 2 Attach the cables to the Preamps.
- 3 Two People must move the Body Hybrid splitter to the rear of the magnet. Maintain maximum distance from magnet at all times.
- 4 Slide the New Hybrid splitter into it saddle.
- 5 Attach the T/R Bias Cable.
- 6 Attach the Cables going to the body coil.
- 7 Attach the Heliax Cables.

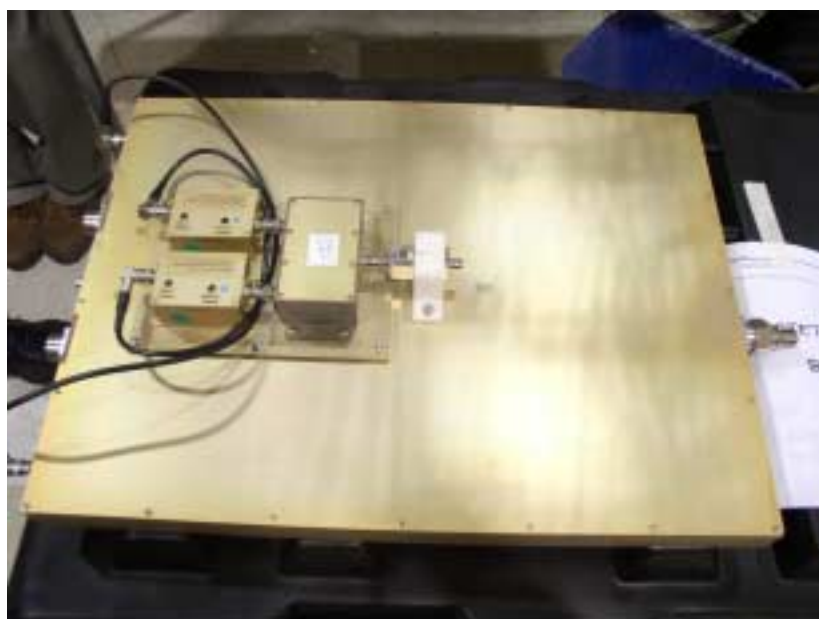

BODY HYBRID/COMBINER ASSEMBLY-READY TO RE-INSTALL
ILLUSTRATION 2-3

2- FINALIZATION

1. Replace All Covers Removed.
2. Complete a Body and Head scan to makes sure the system is functional.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
A	Oct 13, 2003	D. Hofstetter	Initial version