

GE Healthcare

BrightSpeed Select Series Installation Manual

Book 1 of 2: Mechanical Installation

OPERATING DOCUMENTATION



5181248-100
Rev 7
© 2007-2008, General Electric Company
All Right Reserved

Book 1 of 2: Mechanical Installation

Pages 1 - 174

Effectivity:

The information in this manual applies to the following BrightSpeed Select Series CT Systems:

- BrightSpeed Elite Select
- BrightSpeed Edge Select
- BrightSpeed Excel Select

IMPORTANT PRECAUTIONS

LANGUAGE

ПРЕДУПРЕЖДЕНИЕ (BG)	• Това упътване за работа е налично само на английски език. • Ако доставчикът на услугата на клиента изиска друг език, задължение на клиента е да осигури превод. • Не използвайте оборудването, преди да сте се консултирали и разбрали упътването за работа. • Неспазването на това предупреждение може да доведе до нараняване на доставчика на услугата, оператора или пациент в резултат на токов удар или механична или друга опасност.
警告 (ZH-CN)	• 本维修手册仅提供英文版本。 • 如果维修服务提供商需要非英文版本，客户需自行提供翻译服务。 • 未详细阅读和完全理解本维修手册之前，不得进行维修。 • 忽略本警告可能对维修人员，操作员或患者造成触电、机械伤害或其他形式的伤害。
VÝSTRAHA (CS)	• Tento provozní návod existuje pouze v anglickém jazyce. • V případě, že externí služba zákazníkům potřebuje návod v jiném jazyce, je zajištění překladu do odpovídajícího jazyka úkolem zákazníka. • Nesnažte se o údržbu tohoto zařízení, aniž byste si přečetli tento provozní návod a pochopili jeho obsah. • V případě nedodržování této výstrahy může dojít k poranění pracovníka prodejního servisu, obslužného personálu nebo pacientů vlivem elektrického proudu, respektive vlivem mechanických či jiných rizik.
ADVARSEL (DA)	• Denne servicemanual findes kun på engelsk. • Hvis en kundes tekniker har brug for et andet sprog end engelsk, er det kundens ansvar at sørge for oversættelse. • Forsøg ikke at servicere udstyret medmindre denne servicemanual har været konsulteret og er forstået. • Manglende overholdelse af denne advarsel kan medføre skade på grund af elektrisk, mekanisk eller anden fare for teknikeren, operatøren eller patienten.
WAARSCHUWING (NL)	• Deze onderhoudshandleiding is enkel in het Engels verkrijgbaar. • Als het onderhoudspersoneel een andere taal vereist, dan is de klant verantwoordelijk voor de vertaling ervan. • Probeer de apparatuur niet te onderhouden voordat deze onderhoudshandleiding werd geraadpleegd en begrepen is. • Indien deze waarschuwing niet wordt opgevolgd, zou het onderhoudspersoneel, de operator of een patiënt gewond kunnen raken als gevolg van een elektrische schok, mechanische of andere gevaren.

WARNING (EN)	• This Service Manual is available in English only. • If a customer's service provider requires a language other than English, it is the customer's responsibility to provide translation services. • Do not attempt to service the equipment unless this service manual has been consulted and is understood. • Failure to heed this warning may result in injury to the service provider, operator, or patient, from electric shock or from mechanical or other hazards.
HOIATUS (ET)	• Käesolev teenindusjuhend on saadaval ainult inglise keeles. • Kui klienditeeninduse osutaja nõuab juhendit inglise keelest erinevas keeles, vastutab klient tõlketeenuse osutamise eest. • Ärge üritage seadmeid teenindada enne eelnevalt käesoleva teenindusjuhendiga tutvumist ja sellest aru saamist. • Käesoleva hoiatuse eiramine võib põhjustada teenuseosutaja, operaatori või patsiendi vigastamist elektrilöögi, mehaanilise või muu ohu tagajärjel.
VAROITUS (FI)	• Tämä huolto-ohje on saatavilla vain englanniksi. • Jos asiakkaan huoltohenkilöstö vaatii muuta kuin englanninkielistä materiaalia, tarvittavan käännöksen hankkiminen on asiakkaan vastuulla. • Älä yritä korjata laitteistoa ennen kuin olet varmasti lukenut ja ymmärtänyt tämän huolto-ohjeen. • Mikäli tätä varoitusta ei noudateta, seurauksena voi olla huoltohenkilöstön, laitteiston käyttäjän tai potilaan vahingoittuminen sähköiskun, mekaanisen vian tai muun vaaratilanteen vuoksi.
ATTENTION (FR)	• Ce manuel de service n'est disponible qu'en anglais. • Si le technicien du client a besoin de ce manuel dans une autre langue que l'anglais, c'est au client qu'il incombe de le faire traduire. • Ne pas tenter d'intervenir sur les équipements tant que le manuel service n'a pas été consulté et compris • Le non-respect de cet avertissement peut entraîner chez le technicien, l'opérateur ou le patient des blessures dues à des dangers électriques, mécaniques ou autres.
WARNUNG (DE)	• Diese Serviceanleitung existiert nur in Englischer Sprache. • Falls ein fremder Kundendienst eine andere Sprache benötigt, ist es aufgabe des Kunden für eine Entsprechende Übersetzung zu sorgen. • Versuchen Sie nicht diese Anlage zu warten, ohne diese Serviceanleitung gelesen und verstanden zu haben. • Wird diese Warnung nicht beachtet, so kann es zu Verletzungen des Kundendiensttechnikers, des Bedieners oder des Patienten durch Stromschläge, Mechanische oder Sonstige gefahren kommen.

ΠΡΟΕΙΔΟΠΟΙΗΣΗ (EL)	• Το παρόν εγχειρίδιο σέρβις διατίθεται στα αγγλικά μόνο. • Εάν το άτομο παροχής σέρβις ενός πελάτη απαιτεί το παρόν εγχειρίδιο σε γλώσσα εκτός των αγγλικών, αποτελεί ευθύνη του πελάτη να παρέχει υπηρεσίες μετάφρασης. • Μην επιχειρήσετε την εκτέλεση εργασιών σέρβις στον εξοπλισμό εκτός εάν έχετε συμβουλευτεί και έχετε κατανοήσει το παρόν εγχειρίδιο σέρβις. • Εάν δε λάβετε υπόψη την προειδοποίηση αυτή, ενδέχεται να προκληθεί τραυματισμός στο άτομο παροχής σέρβις, στο χειριστή ή στον ασθενή από ηλεκτροπληξία, μηχανικούς ή άλλους κινδύνους.
FIGYELMEZTETÉS (HU)	• Ezen karbantartási kézikönyv kizárólag angol nyelven érhető el. • Ha a vevő szolgáltatója angoltól eltérő nyelvre tart igényt, akkor a vevő felelőssége a fordítás elkészítése. • Ne próbálja elkezdni használni a berendezést, amíg a karbantartási kézikönyvben leírtakat nem értelmezték. • Ezen figyelmeztetés figyelmen kívül hagyása a szolgáltató, működtető vagy a beteg áramütés, mechanikai vagy egyéb veszélyhelyzet miatti sérülését eredményezheti.
AÐVÖRUN (IS)	• Þessi þjónustuhandbók er eingöngu áanleg á ensku. • Ef að þjónustuveitandi viðskiptamanns þarfnast annas tungumáls en ensku, er það skylda viðskiptamanns að skaffa tungumálapjónustu. • Reynið ekki að afgreiða tækið nema að þessi þjónustuhandbók hefur verið skoðuð og skilin. • Brot á sinna þessari aðvörun getur leitt til meiðsla á þjónustuveitanda, stjórnaða eða sjúklings frá raflosti, vélrænu eða öðrum áhættum.
AVVERTENZA (IT)	• Il presente manuale di manutenzione è disponibile soltanto in inglese. • Se un addetto alla manutenzione richiede il manuale in una lingua diversa, il cliente è tenuto a provvedere direttamente alla traduzione. • Si proceda alla manutenzione dell'apparecchiatura solo dopo aver consultato il presente manuale ed averne compreso il contenuto. • Il non rispetto della presente avvertenza potrebbe far compiere operazioni da cui derivino lesioni all'addetto, alla manutenzione, all'utilizzatore ed al paziente per folgorazione elettrica, per urti meccanici od altri rischi.
警告 (JA)	• このサービスマニュアルには英語版しかありません。 • サービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。 • このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないでください。 • この警告に従わない場合、サービスを担当される方、操作員あるいは患者さんが、感電や機械的又はその他の危険により負傷する可能性があります。
경고 (KO)	• 본 서비스 지침서는 영어로만 이용하실 수 있습니다 . • 고객의 서비스 제공자가 영어 이외의 언어를 요구할 경우 , 번역 서비스를 제공하는 것은 고객의 책임입니다 . • 본 서비스 지침서를 참고했고 이해하지 않는 한은 해당 장비를 수리하려고 시도하지 마십시오 . • 이 경고에 유의하지 않으면 전기 쇼크 , 기계상의 혹은 다른 위험으로부터 서비스 제공자 , 운영자 혹은 환자에게 위해를 가할 수 있습니다 .

BRDINJUMS (LV)	• Šī apkalpes rokasgrāmata ir pieejama tikai angļu valodā. • Ja klienta apkalpes sniedzējam nepieciešama informācija citā valodā, nevis angļu, klienta pienākums ir nodrošināt tulkošanu. • Neveiciet aprīkojuma apkalpi bez apkalpes rokasgrāmatas izlasīšanas un saprašanas. • Šī brīdinājuma neievērošana var radīt elektriskās strāvas trieciena, mehānisku vai citu risku izraisītu traumu apkalpes sniedzējam, operatoram vai pacientam.
ĮSPĖJIMAS (LT)	• Šis eksploatavimo vadovas yra prieinamas tik anglų kalba. • Jei kliento paslaugų tiekėjas reikalauja vadovo kita kalba – ne anglų, numatyti vertimo paslaugas yra kliento atsakomybė. • Nemėginkite atlikti įrangos techninės priežiūros, nebent atsižvelgėte į šį eksploatavimo vadovą ir jį supratote. • Jei neatkreipsite dėmesio į šį perspėjimą, galimi sužalojimai dėl elektros šoko, mechaninių ar kitų pavojų paslaugų tiekėjui, operatoriui ar pacientui.
ADVARSEL (NO)	• Denne servicehåndboken finnes bare på engelsk. • Hvis kundens serviceleverandør trenger et annet språk, er det kundens ansvar å sørge for oversettelse. • Ikke forsøk å reparere utstyret uten at denne servicehåndboken er lest og forstått. • Manglende hensyn til denne advarselen kan føre til at serviceleverandøren, operatøren eller pasienten skades på grunn av elektrisk støt, mekaniske eller andre farer.
OSTRZEŻENIE (PL)	• Niniejszy podręcznik serwisowy dostępny jest jedynie w języku angielskim. • Jeśli dostawca usług klienta wymaga języka innego niż angielski, zapewnienie usługi tłumaczenia jest obowiązkiem klienta. • Nie próbować serwisować wyposażenia bez zapoznania się i zrozumienia niniejszego podręcznika serwisowego. • Niezastosowanie się do tego ostrzeżenia może spowodować urazy dostawcy usług, operatora lub pacjenta w wyniku porażenia elektrycznego, zagrożenia mechanicznego bądź innego.
ATENÇÃO (PT-BR)	• Este manual de assistência técnica só se encontra disponível em inglês. • Se qualquer outro serviço de assistência técnica, que não a gems, solicitar estes manuais noutro idioma, é da responsabilidade do cliente fornecer os serviços de tradução. • Não tente reparar o equipamento sem ter consultado e compreendido este manual de assistência técnica. • O não cumprimento deste aviso pode por em perigo a segurança do técnico, operador ou paciente devido a choques elétricos, mecânicos ou outros.

ATENÇÃO (PT-PT)	• Este manual de assistência técnica só se encontra disponível em inglês. • Se qualquer outro serviço de assistência técnica, que não a gems, solicitar estes manuais noutra linguagem, é da responsabilidade do cliente fornecer os serviços de tradução. • Não tente reparar o equipamento sem ter consultado e compreendido este manual de assistência técnica. • O não cumprimento deste aviso pode colocar em perigo a segurança do técnico, do operador ou do paciente devido a choques eléctricos, mecânicos ou outros.
ATENȚIE (RO)	• Acest manual de service este disponibil numai în limba engleză. • Dacă un furnizor de servicii pentru clienți necesită o altă limbă decât cea engleză, este de datoria clientului să furnizeze o traducere. • Nu încercați să reparați echipamentul decât ulterior consultării și înțelegerii acestui manual de service. • Ignorarea acestui avertisment ar putea duce la rănirea depănatorului, operatorului sau pacientului în urma pericolelor de electrocutare, mecanice sau de altă natură.
ОСТОРОЖНО! (RU)	• Данное руководство по обслуживанию предлагается только на английском языке. • Если сервисному персоналу клиента необходимо руководство не на английском, а на каком-то другом языке, клиенту следует самостоятельно обеспечить перевод. • Перед обслуживанием оборудования обязательно обратитесь к данному руководству и поймите изложенные в нем сведения. • Несоблюдение требований данного предупреждения может привести к тому, что специалист по обслуживанию, оператор или пациент получат удар электрическим током, механическую травму или другое повреждение.
UPOZORENJE (SR)	• Ovo servisno uputstvo je dostupno samo na engleskom jeziku. • Ako klijentov serviser zahteva neki drugi jezik, klijent je dužan da obezbedi prevodilacke usluge. • Ne pokušavajte da opravite uređaj ako niste pročitali i razumeli ovo servisno uputstvo. • Zanemarivanje ovog upozorenja može dovesti do povređivanja serviser, rukovaoca ili pacijenta usled strujnog udara ili mehanickih i drugih opasnosti.
UPOZORNENIE (SK)	• Tento návod na obsluhu je k dispozícii len v angličtine. • Ak zákazníkovi poskytovateľ služieb vyžaduje iný jazyk ako angličtinu, poskytnutie prekladateľských služieb je zodpovednosťou zákazníka. • Nepokúšajte sa o obsluhu zariadenia skôr, ako si neprečítate návod na obsluhu a neporozumiete mu. • Zanedbanie tohto upozornenia môže vyústiť do zranenia poskytovateľa služieb, obsluhujúcej osoby alebo pacienta elektrickým prúdom, do mechanického alebo iného nebezpečenstva.

ATENCION (ES)	• Este manual de servicio sólo existe en inglés. • Si el encargado de mantenimiento de un cliente necesita un idioma que no sea el inglés, el cliente deberá encargarse de la traducción del manual. • No se deberá dar servicio técnico al equipo, sin haber consultado y comprendido este manual de servicio. • La no observancia del presente aviso puede dar lugar a que el proveedor de servicios, el operador o el paciente sufran lesiones provocadas por causas eléctricas, mecánicas o de otra naturaleza.
VARNING (SV)	• Den här servicehandboken finns bara tillgänglig på engelska. • Om en kunds servicetekniker har behov av ett annat språk än engelska ansvarar kunden för att tillhandahålla översättningstjänster. • Försök inte utföra service på utrustningen om du inte har läst och förstår den här servicehandboken. • Om du inte tar hänsyn till den här varningen kan det resultera i skador på serviceteknikern, operatören eller patienten till följd av elektriska stötar, mekaniska faror eller andra faror.
DIKKAT (TR)	• Bu servis kilavuzunun sadece ingilizcesi mevcuttur. • Eğer müşteri teknisyeni bu kilavuzu ingilizce dışında bir başka lisandan talep ederse, bunu tercüme ettirmek müşteriye düşer. • Servis kilavuzunu okuyup anlamadan ekipmanlara müdahale etmeyiniz. • Bu uyarıya uyulmaması, elektrik, mekanik veya diğer tehlikelerden dolayı teknisyen, operatör veya hastanın yaralanmasına yol açabilir.

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation "damage in shipment" written on all copies of the freight or express bill before delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage MUST be reported to the carrier immediately upon discovery, or in any event, within 14 days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this 14 day period.

To file a report:

- Call 1-800-548-3366 and use option 8.
- Fill out a report on <http://egems.med.ge.com/edq/home.jsp>
- Contact your local service coordinator for more information on this process.

Rev. June 13, 2006

CERTIFIED ELECTRICAL CONTRACTOR STATEMENT

All electrical Installations that are preliminary to positioning of the equipment at the site prepared for the equipment shall be performed by licensed electrical contractors. In addition, electrical feeds into the Power Distribution Unit shall be performed by licensed electrical contractors. Other connections between pieces of electrical equipment, calibrations and testing shall be performed by qualified GE Medical personnel. The products involved (and the accompanying electrical installations) are highly sophisticated, and special engineering competence is required. In performing all electrical work on these products, GE will use its own specially trained field engineers. All of GE's electrical work on these products will comply with the requirements of the applicable electrical codes.

The purchaser of GE equipment shall only utilize qualified personnel (i.e., GE's field engineers, personnel of third-party service companies with equivalent training, or licensed electricians) to perform electrical servicing on the equipment.

IMPORTANT...X-RAY PROTECTION

X-ray equipment if not properly used may cause injury. Accordingly, the instructions herein contained should be thoroughly read and understood by everyone who will use the equipment before you attempt to place this equipment in operation. The General Electric Company, Medical Systems Group, will be glad to assist and cooperate in placing this equipment in use.

Although this apparatus incorporates a high degree of protection against x-radiation other than the useful beam, no practical design of equipment can provide complete protection. Nor can any practical design compel the operator to take adequate precautions to prevent the possibility of any persons carelessly exposing themselves or others to radiation.

It is important that anyone having anything to do with x-radiation be properly trained and fully acquainted with the recommendations of the National Council on Radiation Protection and Measurements as published in NCRP Reports available from NCRP Publications, 7910 Woodmont Avenue, Room 1016, Bethesda, Maryland 20814, and of the International Commission on Radiation Protection, and take adequate steps to protect against injury.

The equipment is sold with the understanding that the General Electric Company, Medical Systems Group, its agents, and representatives have no responsibility for injury or damage which may result from improper use of the equipment.

Various protective materials and devices are available. It is urged that such materials or devices be used.

LITHIUM BATTERY CAUTIONARY STATEMENTS



CAUTION
Risk of
Explosion

Danger of explosion if battery is incorrectly replaced. Replace only with the same or equivalent type recommended by the manufacturer. Discard used batteries according to the manufacturer's instructions.



ATTENTION
Danger
d'Explosion

Il y a danger d'explosion s'il y a remplacement incorrect de la batterie. Remplacer uniquement avec une batterie du même type ou d'un type recommandé par le constructeur. Mettre au rebut les batteries usagées conformément aux instructions du fabricant.

OMISSIONS & ERRORS

Customers, please contact your GE Sales or Service representatives.

GE personnel, please use the GEMS CQA Process to report all omissions, errors, and defects in this publication.

Table of Contents: Book 1

Preface

Publication Conventions 19

Section 1.0

Safety & Hazard Information 19

1.1 Text and Character Representation..... 19

1.2 Graphical Representation 20

Section 2.0

Publication Conventions 21

2.1 General Paragraph and Character Styles..... 21

2.2 Page Layout..... 21

2.3 Computer Screen Output/Input Character Styles 22

2.4 Buttons, Switches and Keyboard Inputs (Hard & Soft Keys)..... 22

Chapter 1

Position Subsystems 23

Section 1.0

Installer/FE Notices 23

1.1 Shipping, Warehouse, and Transportation Warning..... 23

1.2 International Shipments 23

1.3 On Site Warning (Mechanical Installation Team) 23

Section 2.0

Install Table/Gantry Introduction 23

2.1 Skill Set Required 24

2.2 Floor and Room Preparation 24

2.2.1 Preparation 24

2.2.2 Flooring..... 25

2.3 Overview..... 25

2.4 Pre-Installation Template..... 25

2.5 Required Common Tools and Supplies 26

Section 3.0

Delivery Procedure..... 28

3.1 System Transportation - Temperature Extremes..... 28

3.2 Stored Systems 28

3.3 Working with the Mover 28

3.4 Floor Protection 28

3.5 Equipment Delivery Route 28

3.6 Removing Gantry Dollies and Covers..... 28

3.7 Damage In Transportation 29

3.8 A1 Breaker..... 29

3.9 Installation Support Kits 29

3.10 Installation Conditions..... 29

Section 4.0

Install and Level Table/Gantry 30

4.1	Establish the Room Layout	30
4.2	Position the Gantry.....	31
4.2.1	Gantry Prep - For Access Greater Than 28"	31
4.2.2	Gantry Prep - For Access Less Than 28"	31
4.2.3	Gantry Positioning - All Sites	32
4.3	Level the Gantry	34
4.4	Gantry Bearing Gap Inspection	36
4.4.1	Personnel Requirements.....	36
4.4.2	Tools and Test Equipment	36
4.4.3	Damage Indicators	36
4.4.4	Porcedure.....	37
4.4.5	Finalization - Mechanical Installers	38
4.4.6	FE Service Action Required	39
4.4.7	FE Inspection Completion	39
4.5	Table Prep and Setup	40
4.6	Cradle Center Procedure	42
4.7	Center line on Cradle	43
4.8	Install Gantry Alignment Laser and Bracket.....	44
4.9	Install the Table Cradle Laser Alignment Plates	47
4.10	Level and Center the Table to the Gantry	48
4.10.1	Required Tools.....	48
4.10.2	Alignment Conditions	48
4.10.3	Alignment Specification	48
4.10.4	Procedures.....	49
4.11	Drill the Anchor Holes	50
4.11.1	Notes to Mechanical Installers	50
4.11.1.1	Note 1: Alternate Anchoring	50
4.11.1.2	Note 2: Non-Concrete Floors	50
4.11.1.3	Note 3: GE Notification	50
4.11.2	Gantry Drilling	51
4.11.3	Table Drilling and Anchoring	52
4.12	Gantry & Table Redundant Anchor Holes	55
4.13	Gantry Anchoring	56
4.14	Alignment Recheck	57
Section 5.0		
Rear Entry Cable Box.....		58
Section 6.0		
Install Table Footswitch Assembly		58
Section 7.0		
Remove Gantry Tilt Bracket.....		61
Section 8.0		
Position the Power Distribution Unit		62
Section 9.0		
Install Operator Console Without FWS Table		64
9.1	Unpack Console	64
9.2	Remove Console Covers	66
9.3	Adjust Table Top Height and Position Console.....	66
9.4	Attach Keyboard Table Top	67

9.5	Monitor and SCSI Tower Placement	68
-----	--	----

Section 10.0

Install Operator Console With FWS Table	69
--	-----------

10.1	Unpack Console	69
10.2	Remove Console Covers.....	70
10.3	Replace Console Table Top With Black Table Top	70
10.4	Install FWS table and Monitor Arms	71
10.5	Peripherals Placement.....	71

Section 11.0

Seismic Mounting.....	72
------------------------------	-----------

11.1	Console.....	72
11.2	Power Distribution Unit	72

Chapter 2

Power, Ground & Interconnect Cables	73
--	-----------

Section 1.0

Introduction	73
---------------------------	-----------

1.1	System Component Identification	74
1.2	Cable Color Identifiers	74

Section 2.0

System Interconnect Diagram.....	76
---	-----------

Section 3.0

Contractor Connections	78
-------------------------------------	-----------

Section 4.0

Console Connections	79
----------------------------------	-----------

4.1	SCIM, Keyboard, Trackball (Option) & Mouse Installation	79
4.2	Connecting the SCSI Tower	81
4.3	Connecting the Monitor.....	81
4.4	Power Panel Connections	82
4.5	Installing a DASM (Hardware Option).....	84
4.6	Modem Option	85
4.7	LAN Connections.....	85
4.8	Console Interconnect.....	86

Section 5.0

Install Options	88
------------------------------	-----------

5.1	Install Optional Bar Code Reader	88
5.2	Install Optional Remote Monitor	88
5.3	Install Injector Option	88
5.4	Install The IVY Monitor and Stand	88
5.5	Customer Accessories (Head Holders and Extender)	88
5.6	Install Service Cabinet.....	88
5.7	Install UPS	88

Section 6.0

Gantry Cable Connections	90
---------------------------------------	-----------

Section 7.0	
Table Connections.....	93
Section 8.0	
PDU Cable Connections & Configuration	94
8.1 NGPDU	94
8.1.1 Panel - 380 - 480VAC Mains “TS1” Input Power Connection	95
8.1.2 Panel - Circuit Breakers	96
8.1.3 Transformer (480VAC) Taps	97
8.1.4 HVDC Connection.....	97
8.1.5 440V Connection.....	98
8.1.6 Gantry & Console Power Connections (120V).....	98
8.1.7 Console Power Cable Re-termination.....	98
8.1.8 PDU Control Cable.....	99
8.1.9 System Ground Connection	99
8.1.10 Warning Light & Door Interlock Connections	99
8.1.10.1 Warning Light Configuration & Connection.....	99
8.1.10.2 Door Interlock Connections.....	100
Section 9.0	
System Ground Connections	102
 Chapter 3	
System Continuity & Ground Checks	105
Section 1.0	
System Continuity & Ground Checks	
(Mechanical Contractor Task)	105
1.1 Tools Required.....	105
1.2 Procedure.....	105
Section 2.0	
Site Ground Continuity Check.....	108
Section 3.0	
Axial Head Holder Shim Installation	109
3.1 Overview	109
3.2 Requirements.....	109
3.3 Required Conditions.....	110
3.4 Procedure.....	111
3.5 Finalization	111
Section 4.0	
Mechanical Installation Completion Section.....	112
 Appendix A	
Removal & Installation of Covers	113
Section 1.0	
Gantry Side Covers	113
1.1 Side Cover Removal	113
1.2 Side Cover Installation	114

Section 2.0	
Gantry Top Covers	114
2.1 Top Cover Removal	114
2.2 Top Cover Installation	115
Section 3.0	
Gantry Front Cover	116
3.1 Front Cover Dolly Setup	116
3.2 Removal.....	117
3.3 Installation.....	122
Section 4.0	
Gantry Rear Cover.....	124
4.1 Removal.....	124
4.2 Installation.....	125
Section 5.0	
Gantry Cover Alignment Adjustment Guide	126
5.1 Overview	126
5.2 Tools Required	126
5.3 Adjustment Guide	126
5.3.1 Front/Rear Cover adjustment	126
5.3.2 Scan Window	127
5.3.3 Top Cover Adjustment	127
5.3.4 Side Covers	128
5.3.5 Finish	128
Section 6.0	
Gantry Scan Window	128
6.1 Remove Scan Window	128
6.2 Install Scan Window	128
Section 7.0	
Gantry Base Covers	129
Section 8.0	
Table Covers.....	131
8.1 Table Base Cover	131
8.2 Flapper Installation	131
8.3 Skirt Cover installation	132
8.4 Bottom cover installation.....	133

Appendix B

FWS Assembly and Adjustment **135**

Section 1.0	
Introduction	135
Section 2.0	
Procedure.....	136
2.1 Assemble Worksurface (For Part 5160735)	136
2.2 Assemble Worksurface (For Part 5168666)	139
2.3 Install Monitors.....	145

2.4	Monitor Adjustment	150
2.5	Install Seismic Kit (For Part 5160735).....	151
2.6	Install Seismic Kit (For Part 5168666).....	152

Appendix C

Pictorial Representation of Required Tools 153

Appendix D

Option Cables 157

Appendix E

Regulatory Clearance Quick Reference Guide..... 159

Section 1.0	
	Regulatory Code Description 159
Section 2.0	
	Terms and Definitions 160
Section 3.0	
	Regulatory Minimum Working Clearance by Major Subsystem..... 161
Section 4.0	
	Additional Regulatory Clearance Information 165
4.1	Regulatory Caution 165
4.2	Operational Caution 165
4.3	Recommended Room Size 165
4.4	Typical Room Size 165
4.5	Minimum Room Size 166
Section 5.0	
	System Specifications (BrightSpeed Select) 166
5.1	Rooms w/ Less Than 711mm (28 in) Egress Clearance at Table Foot End 166
Section 6.0	
	How to Measure 167
Section 7.0	
	Minimum Room Size & Requirement Layouts 168
7.1	Typical and Recommended Room Size & Requirement Layouts 170

Revision History

Revision	Date	Reason for change
7	11/04/08	Chapter 1: Section 4.0 updated Table/Gantry Alignment Procedure. Chapter 2: Section 2.0 added system interconnected diagram with AC Outlet Box Section 4.4 added figure of power connection with AC Outlet Box updated Section 4.6 Modem Option for AIO console updated Secion 4.8, added AIO Interconnect updated Section 5.5 for Install Injector Option Section 9.0 updated figure of CT System Ground Connections Chapter 5: updated Section 2.10, added AIO Console Boot-up Flow Chart Section 3.0 corrected Test Item 3.10.2.2-6 and 3.10.2.2-8 values Section 3.4 updated CT System X-Ray ON Indicators, Cautions & Warning Labels Section 3.11.3 Deleted Alignment Light Characterization, it is duplicated with 3.9.1. Chapter 6: Section 3.0 updated the procedure with Laser Alignment Tool. Chapter 7: Section 2.0 Deleted install Injector and UPS, they are duplicated with Chapter 2, Section 5. added Section 11 for Modem setup
6	06/ 01/08	Chapter 2: updated Section 4.6 Modem Option for AIO console updated Secion 4.8, added AIO Interconnect added Section 5.0 Install Options Appendix B: added new FWS (part# 5168666) installation procedure. Chapter 5: updated Section 2.10, added AIO Console Boot-up Flow Chart
5	04/21/08	Chapter 2: added Section 5.0 Install Options Appendix E: updated room size descriptions Chapter 5: updated Section 3.7 Verify Service Cabinet Installation updated Section 3.9, 3.10 for Table/Gantry interference test specifications. Customer Configuration Preference for Language requirements: Updated Chapter 5, Section 2.6.2, Step 7c Procedure, and Section 3.5 Warning Labels. Also updated GE form e4879 to record information.

Revision	Date	Reason for change
4	10/15/07	Important Precautions: Added Bulgarian warnings and damage in transportation Chapter 1: added Section 1.0 installer/FE notices updated Section 2.1 Training Process updated Section 2.5 common tools updated Section 4.3 Table Prep and Setup updated Section 4.5 procedure, added Cable entrance width Chapter 2: updated Section 4.4 injector connection port updated Section 4.7 connection cable Chapter 3: updated Section 4.0 Mechanical Installation Completion Section Appendix B: added Section 2.4 install Seismic kit Appendix E: removed Gantry bearing gap inspection procedures; add related procedures to Chapter 1 section 4.4 Chapter 4: updated Section 4.3 phantom added Section 6.0 FE Workflow added Section 7.0 Checklist for Completed Installation added Section 8.0 GE and Regulatory Forms Chapter 5: updated Section 2.3 system state DVD location updated Section 2.7 Check/Set Date and Time Chapter 6: updated specs in Section 6.4 Chapter 7: updated Section 2.0 UPS installation Chapter 8: updated Patient Touch Leakage Test. added the CT System Chassis Leakage Test Chapter 9: removed this chapter ("Forms" chapter). Forms now reside on the Service CD. Appendix F: updated Regulatory Clearance Quick Reference Guide
3	6/1/07	Chapter 1: updated Gantry/Table alignment tools installation procedure updated Table anchoring procedure Chapter 2: updated Trackball as option Chapter 9: updated Table/Gantry Alignment procedure Checklist
2	4/17/07	Chapter 1: updated gantry drill bit. added Installation Conditions. added Limited Access Installation instructions. updated redundant anchor locations. Appendix E: added Gantry bearing gap inspection procedures. Appendix F: added Regulatory Clearance Quick Reference Guide.
1	4/11/07	Initial Release.

Preface

Publication Conventions

Please become familiar with the conventions used within this publication before proceeding.

Section 1.0 Safety & Hazard Information

1.1 Text and Character Representation

Within this publication, different paragraph and character styles have been used to indicate potential hazards. Paragraph prefixes, such as hazard, caution, danger and warning, are used to identify important safety information. Text (Hazard) styles are applied to the paragraph contents that are applicable to each specific safety statement. Words describe the type of potential hazard that may be encountered and are placed immediately before the paragraph it modifies. Safety information will normally include:

- Type of potential hazard
- Nature of potential injury
- Causative condition
- How to avoid or correct the causative condition

EXAMPLES OF HAZARD STATEMENTS USED

A few examples are provided that have been adapted from GEMS' global document standard (2119696-100). They include paragraph prefixes and modified text styles.



CAUTION
Pinch Points
Loss of Data
Sharp Objects

Caution is used when a hazard exists that can or could cause minor injury to self or others if instructions are ignored. They include for example:

- Loss of critical patient data
- Crush or pinch points
- Sharp objects



DANGER
EXCESSIVE
VOLTAGE
CRUSH
POINT

DANGER IS USED WHEN A HAZARD EXISTS THAT WILL CAUSE SEVERE PERSONAL INJURY OR DEATH IF INSTRUCTIONS ARE IGNORED. THEY CAN INCLUDE:

- ELECTROCUTION
- CRUSHING
- RADIATION



WARNING
ROTATING
EQUIPMENT
BARE WIRES

WARNING IS USED WHEN A HAZARD EXISTS WHICH COULD OR CAN CAUSE SERIOUS PERSONAL INJURY OR DEATH IF INSTRUCTIONS ARE IGNORED. THEY CAN INCLUDE:

- Potential for shock
- Exposed wires
- Failure to Tag and lockout system power could allow for un-command motion.

**NOTICE**
Equipment Damage Possible

Notice is used when a hazard is present that can cause property damage but has absolutely no personal injury risk. They can include:

- Disk drive will crash
- Internal mechanical damage, such as to the x-ray tube
- Coasting the rotor through resonance.

It's important that the reader not ignore hazard statements in this document.

1.2 Graphical Representation

Important information will always be preceded by the exclamation point  contained within a triangle, as seen throughout this chapter. In addition to text, several different graphical icons (symbols) may be used to make you aware of specific types of hazards that could possibly cause harm.

ELECTRICAL



LASER



MECHANICAL



HEAT



RADIATION



PINCH



Some others make you aware of specific procedures that should be followed.

AVOID STATIC ELECTRICITY



TAG AND LOCK OUT



WEAR EYE PROTECTION



Section 2.0 Publication Conventions

2.1 General Paragraph and Character Styles

Prefixes are used to highlight important non-safety related information. Paragraph prefixes (such as Purpose, Example, Comment and Note) are used to identify important but non-safety related information. Text styles are also applied to text within each paragraph modified by the specific prefix.

EXAMPLES OF PREFIXES USED FOR GENERAL INFORMATION


Purpose: Introduces and provides meaning as to the information contained within the chapter, section or subsection, such as used at the beginning this chapter for example.

Note: Conveys information that should be considered important to the reader.

Example: Used to make the reader aware that the paragraph(s) that follow are examples of information possibly stated previously.

Comment: Represents "additional" information that may or may not be relevant.

2.2 Page Layout

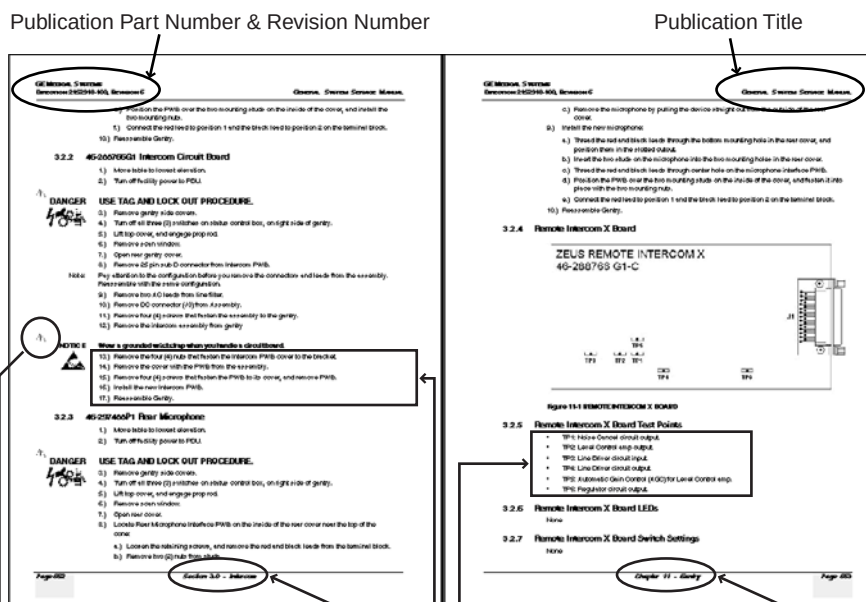

The current section and its title are always shown in the footer of the left (even) page.

An exclamation point in a triangle is used to indicate important information to the user.

Paragraphs preceded by Alphanumeric (e.g. numbers) characters is information that must be followed in a specific order.

The current chapter and its title are always shown in the footer of the right (odd) page.

Paragraphs preceded by alphanumeric (e.g. numbers) characters is information that has no specific order.

Headers and footers in this publication are designed to allow you to quickly identify your location. The document's part number and revision number appears in every header on every page. Odd

numbered page footers indicate the current chapter, its title, and current page number. Even page footers show the current section and its title, as well current page number.

2.3 Computer Screen Output/Input Character Styles

Within this publication different character styles are used to indicate computer input and output text. Character (input, output, and variable) styles are used and applied to the text within a paragraph so as to indicate directions. Computer screen output and input is also formatted using mono (fixed width) spaced fonts.

Example: This paragraph denotes computer screen fixed output. It's output is fixed from the sense that it does not vary from application to application. It's the most commonly used style used to indicate filenames, paths, and text.

Example: *This paragraph denotes computer screen output that is variable. Its output varies from application to application. Variable output is sometimes found placed between greater than and lesser than operators. For example: <variable_ouput>*

Example: **This paragraph denotes fixed input. It's typed input that will not vary from application to application. Fixed text the user is required to supply as input.**

Example: ***This paragraph denotes computer input that can vary from application to application. Variable text the user is required to supply as input. Variable input sometimes is placed between greater than and lesser than operators. For example: <variable_input>. In these cases, the (<>) operators are dropped prior to input. Exceptions are noted in the text.***

2.4 Buttons, Switches and Keyboard Inputs (Hard & Soft Keys)

Different character styles are used to indicate actions requiring the reader to press either a hard or soft button, switch, or key. Physical hardware, such as buttons and switches, are called hard keys because they are hard wired or mechanical in nature. A keyboard or on/off switch would be a hard key. Software or computer generated buttons are called soft keys because they are software generated. Software driven menu buttons are an example of such keys. Soft and hard keys are represented differently in this publication.

Example: A power switch **ON/OFF** or a keyboard key like **ENTER** is indicated by applying a character style that uses both over and under-lined bold text that is bold. This is a hard key.

Example: Whereas the computer **MENU** button that you would click with your mouse or touch with your hand uses over and under-lined regular text. This is a soft key.

Chapter 1

Position Subsystems



NOTICE

Potential for Data Loss and/or Equipment Damage.

To prevent potential data loss and equipment damage, please do the following:

- When instructed, record data collected from the procedures in this chapter on GE Form e4879. For more information about this form, see [Section 8.0](#) of Chapter 4.
- Only use the Installation manual that arrives with your system. Any other revisions of this manual may not exactly match your system.

Section 1.0

Installer/FE Notices

1.1 Shipping, Warehouse, and Transportation Warning

- This gantry is designed to be moved using the shipping dollies and should not be lifted or moved using a lift truck.
- Do not hoist gantry or table using dollies.

1.2 International Shipments

- Dollies must be used to remove the gantry from the shipping skid and to transport the gantry to the customer's site.
- If lifting is required, refer to the Pre-Installation Manual for instructions.

1.3 On Site Warning (Mechanical Installation Team)

This system requires a gantry bearing gap inspection before electrical calibration is started. See [Gantry Bearing Gap Inspection, on page 36](#).

Section 2.0

Install Table/Gantry Introduction

This chapter describes how to mount, position, and level the CT Scanner subsystems.

Note: Before you start the installation, make sure the site preparation complies with conditions and instructions found in the pre-installation manual. Failure to comply will result in excessive installation delay and potential increased, unrecoverable installation costs. This product is designed to meet specific mechanical installation standards, which should be reviewed prior to installing this system.

2.1 Skill Set Required

NOTICE All Mechanical Installers must meet these minimum requirements.



MECHANICAL

Installation	Safety	OJT with Trainer
Must be trained as FE or Contractor.	Complete all training courses.	Complete training installation with trainer.
Must complete all training Courses.	My Learning verification.	
My Learning verification.		
Tools	Access	Communication
Must have all required tools as listed in this manual.	Must have access to some GE Healthcare web sites and tools.	Must be able to complete all required paperwork and complete iTracs as required.

ELECTRICAL

In-Residence Training.	Safety	OJT with Trainer
Must be trained as FE or Contractor.	Complete all training courses.	Complete training installation with trainer.
My Learning verification.	My Learning verification.	

Only people experienced and trained in the installation process of a GE CT Scanner should perform the work described in this document. A person doing the installation should have:

- Installed three (3) Lightspeed 5.X series or BrightSpeed (Select) series products as a helper prior to leading an installation.
- A working understanding of the mechanical installation process.
- Completed Environmental Health and Safety (EHS) training, completed Hazcom training, completed Eye & Personal Protection Safety, completed Blood Borne Pathogens training, and must understand and apply lockout/ tagout processes
- Professional communication skills needed to communicate with customer any changes or additional requirements needed to complete the installation project.
- Must be trained to complete site-ready inspections.
- (Installation leader) Must have completed project management, leadership training, and all required training prior to becoming a project leader.
- (Mechanical installation teams) Must have basic measurement skills, system layout skills, and the ability to read and understand blueprints.

2.2 Floor and Room Preparation

2.2.1 Preparation

Consult your local GE Sales and service representative about your specific needs. It is the purchaser's (buyer) responsibility to provide an approved support structure and an approved method of mounting. General Electric is not responsible for any failure of the support structure or method of anchoring.

The BrightSpeed Select system has a total floor load of approximately 6260 lbs (2840 kg). About 5000 lbs (2270 kg), including patient (450 lbs (205 kg)), is concentrated in the table-gantry assembly. Refer to the *Pre-installation Manual* for more information.

2.2.2 Flooring

Do not place the scanner on any resilient flooring. Resilient tile or carpeting may slowly yield over a period of time and disturb the alignment of the table to the gantry. Refer to the floor template to determine locations where resilient flooring material should be removed.

Limitations include:

- No part of the floor surface within the table, gantry, or the two interface areas between table and gantry, should be higher than the support areas for the table and gantry.
- The floor structure must withstand the occupied weight of table and gantry, as well as the individual contact area loading of these components.
- The method and placement of anchors, or through bolts, must not reduce the structural strength of the floor.

If you have to remove the gantry covers in order to move the gantry into the room, refer to [Appendix A](#), for the cover removal procedure. Please read the caution statement on [page 113](#) before removing the gantry covers.

2.3 Overview

Procedures in this chapter provide detailed instructions to position, level, and anchor the gantry and table securely for operation. The BrightSpeed Select system uses adjustable leveling pads to support the gantry and table. The gantry has four (4) primary leveling pads located on the gantry base. The table has four (4) pads used for leveling it.

The process you will be following is:

- 1.) Use the room-layout template to determine the general position of the gantry and table.
- 2.) Move the gantry into position.
- 3.) Level gantry.
- 4.) Use the alignment tool to position the table relative to the gantry.
- 5.) Level the table to the gantry, and anchor the system.

Use the template to position the system; however, use the gantry and table to locate and drill the anchor holes. Drill the anchor holes with the system in place.

2.4 Pre-Installation Template

Always use the room-layout template (in two pieces), during installation (see [Table 1-1](#), for part number). The gantry and table will not be properly aligned if existing holes are used. The template shows the location of the gantry and table anchor holes.

This template is shipped with the system. It may also be ordered as FRU.

Table 1-1 Room Layout Templates

System	Part Number
BrightSpeed Select Elite, Edge, Excel	5193991

2.5 Required Common Tools and Supplies

The following tools and supplies are required for installation of the scanner. Please refer to [Appendix C, on page 153](#), for pictorial descriptions of the tools and supplies.

WRENCHES

- Standard and Metric combination wrench sets
- Standard and Metric Hex Key (Allen wrench) sets
- ½" and 3/8" drive torque wrench: 0-100 N-m (0-100 ft.-lb.) Must be calibrated yearly.

SOCKETS AND EXTENSIONS

- 3/8" and ½" drive ratchet wrenches
- ½" drive 3" & two 6" long extensions
- 3/8" drive 12" long extension
- Standard & Metric 3/8" drive socket sets
- ¾" deep well socket 3/8"
- 1", 1-1/8", 1-¼" & 1-½" sockets for ½" drive
- 3/8" drive universal joint
- Metric hex bit set ¼" or 3/8" drive, including:
 - 14mm hex bit 3/8" or ½" drive (14mm ball hex helpful)
 - 10mm hex bit 3/8" drive

SCREW DRIVERS

- Phillips screwdriver set (small, medium, and large)
- Straight blade screwdriver set (small, medium, and large)

DRILL BITS

- Complete set of standard (U.S.) drill bits
- Metric tap set
- 12 mm concrete drill bit (2203081: HILTI, TE-CX 12/17 or equivalent)
- ½" masonry bit, min. 8" long USA – 18" optional (for rear table hole)
- 3" (76mm) hole saw with 1/4" (6mm) masonry bit (to remove flooring)

POWER TOOLS

- 3/8" or ½" drill, cordless or electric
- Reciprocating Saw (Sawzall or equivalent) and assorted blades.
- Hammer Drill & Bit (8" min, 12" max)
- Sears 17740 Shop vacuum or equivalent, with "HEPA" or dry wall dust filter (Sears part number 17918) or equivalent
- 25' Extension power cords

HAND TOOLS

- Ball-Peen Hammer (1lb or 2lb)
- Tongue & Groove Pliers (large)
- Diagonal Cutting Pliers, Large (to cut 1/0 ground)
- Framing Square (e.g., Empire 16" x 24" aluminum square)
- Diagonal Cutting Pliers, Small
- Large pry bar
- 4', 2' & 9" torpedo levels (see [Table 1-2: Recommended Levels](#))
- Laser level (see [Table 1-2: Recommended Levels](#))

Table 1-2 Recommended Levels

9"	Johnson Magnetic Level, model 7500M*
	Johnson Magnetic Level, model 4500
	Stanley Magnetic Level
2'	Johnson Professional Box Beam Level, model 9624*
	Empire Titan Professional Box Beam Level, model 900 series
4' (nominal)	48" Johnson Professional Box Beam Level*
	42" Stanley Contractor Grade Level
	48" Empire Titan Professional Box Beam Level, model 900 series
	48" Stabila Aluminum Box Beam Level, Kit 24816
Laser**	Sears Laser Level Tool (\$40)

* Preferred levels
** Laser level is required for table cradle alignment. Using this tool reduces the alignment procedure by 1 man-hour. New tables have a cradle center line to be used with the laser level.

ELECTRICAL TOOLS

- DVM
- Continuity tester

PERSONAL SAFETY EQUIPMENT

- Safety shoes*
- Safety glasses*
- Gloves
- LOTO Kit (supplied)*
- Hearing Protection*
- 6' & 8 ' Step ladder

* These PPE items are absolutely required for every installation job, with NO exceptions.

SYSTEM CLEANERS

Purchase Locally:

- Glass Cleaner or equivalent
- Scrubbing Bubbles Bathroom Cleaner or equivalent
- Grease Cutting Cleaner or equivalent

GE TOOLS

- 5148193 (VCT alignment Kit)
- 5180876 (TARGET TABLE ALIGNMENT)
- 5272090 (Complete Laser alignment tool for VCT, Brightspeed Select and RT16, which includes 5148193 and 5180876)

Note: A box labeled Installation Support Kit is shipped with each system. It contains paint, masking tape, cleaners, towels, and other materials needed to install this CT scanner.

Section 3.0

Delivery Procedure

3.1 System Transportation - Temperature Extremes

When transporting the CT system, ensure that the system is not exposed to temperatures or humidity outside the following specifications.

Temperature: -29° to +140° F (-34° to +60° C)

Humidity: 5% to 95%

NOTICE **Component Freezing occurs if CT system is exposed to temperatures below 0° F (-18° C) for a period longer than two days.**

Allow a minimum of 12 hours for the CT system to adjust to ambient room temperature, prior to installation.

3.2 Stored Systems

If your system have been stored for more than three months, you will need to complete a visual inspection, looking for damage due to improper storage. Check for the latest software revisions, options, and component changes. Contact the OLC or CT support central for additional information.

3.3 Working with the Mover

Follow the instructions provided by your installation specialist regarding working with equipment movers. Help direct movers on where to place equipment and which items you need first.

Generally movers should move all equipment into the customer room. Door removal and other site changes to move equipment should be done only as directed by the install specialist.

For component sizes and weights, refer to the *BrightSpeed Select Elite, Edge, Excel Pre-Installation* manual.

3.4 Floor Protection

It is suggested that the movers use floor protection. Most equipment movers can provide floor protection during the equipment delivery. Installers should provide floor protection for the room.

3.5 Equipment Delivery Route

Prior to equipment delivery, review the delivery route with the movers. Refer to the installation specialist for any additional delivery instructions.

3.6 Removing Gantry Dollies and Covers

Refer to [Appendix A](#), for the dolly and cover removal procedures. Please read the caution statement on [page 113](#) before removing the gantry covers.

3.7 Damage In Transportation

Check for the following damage, which may have occurred during shipping:

- Equipment damage: paint scrapes and cover damage.
- Damage to hospital property: floors, door frames and walls.

If damage is found—or items are missing in shipment—notify the appropriate service personnel:

- For item(s) missing in shipment or short shipped, contact the install specialist to have the item(s) shipped.
- Report damaged items:
<http://egems.med.ge.com/edq/home.jsp>.

3.8 A1 Breaker

Lock-out and tag-out the A1 breaker now.

Figure 1-1 Sample A1 Breaker



3.9 Installation Support Kits

An Installation Support Kit is shipped with every system. Locate this box now and open it. All included materials are to be used during the installation process. These items are to be left ON SITE, for future service needs.

3.10 Installation Conditions

- 1.) A Final Site Print is REQUIRED.
- 2.) The room size must match the print. Measure the room size. If it does not match the stated size, and is smaller, then check all regulatory clearances. If any regulatory clearance is less than the minimum, then DO NOT continue. Notify the PMI to set up a site escalation.
- 3.) Do Not start the installation process if the site is under construction:
 - In the Room
 - In the Scan AreaRefer to the Pre-Installation Manual for additional details.
- 4.) A customer Anchoring Plan is required, if there is anything other than a 4" (minimum) concrete floor. GEH employees shall only install the anchors supplied with this system.

Section 4.0

Install and Level Table/Gantry

4.1 Establish the Room Layout



NOTICE The floor must meet the minimum floor specification in order to begin this procedure. Refer to the *Pre-Installation Manual* for details.


Use the GE print (developed for your site) to establish the room layout. Make sure all the operating and service clearances shown on the print are observed. Using the supplied template, locate the anchor holes. Make sure they clear structural interferences in the floor.

Clean the area. Free the mounting surface of any material that may interfere with the positioning and leveling of the system.

- 1.) Lay out the two (2) pieces of the floor template (P/N 5193991).
- 2.) Start with the Gantry template—align per the GE print.

- Note: If the floor has existing anchor penetration, the new anchors must be set at least 4" from the nearest penetration.
- 3.) Place the table template over the top of the Gantry template. Align the scan and table center-lines and secure the templates to the floor. Make sure there are no potential clearance issues.
 - 4.) Now, check the level of the floor (See [Figure 1-2](#)) across the templates.

Figure 1-2 Hole Locations

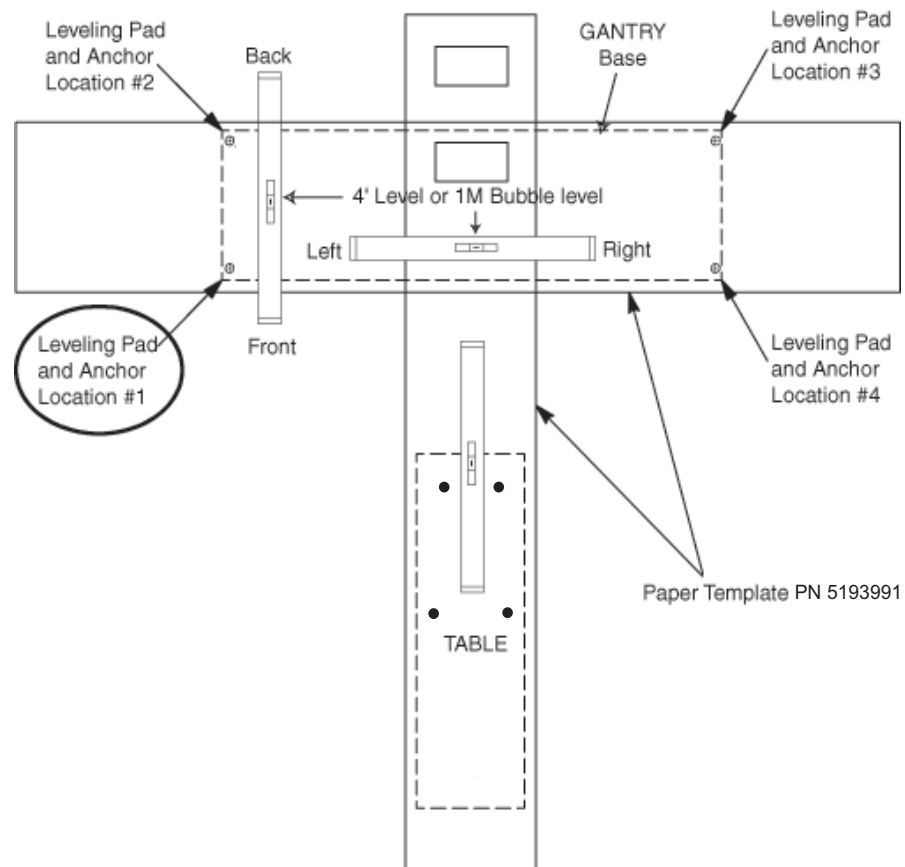

- 5.) Scribe a mark (e.g., use a center punch) at each of the gantry's mounting hole locations (there

are four (4) of these).

- 6.) Using a center punch, mark the four (4) table mounting holes.

NOTICE



Positioning requires cutting 8 holes in the floor.

Before you drill or cut any flooring, make sure the appropriate hospital personnel have approved the location of the table/gantry.

- 7.) Cut tiles (or other resilient flooring) around all holes punched in the template for the gantry and table. Use a utility knife with a heat gun, a 3" hole saw with a 1/4" masonry bit or other adequate tools to cut the flooring.
- 8.) Some sites may require sealing of the floor penetrations after the flooring is removed.

- Note:
- Use RTV or other sealant to seal the floor covering, as necessary.
 - All documentation in this manual is based on mounting the table / gantry on a concrete floor only.

4.2 Position the Gantry

4.2.1 Gantry Prep - For Access Greater Than 28"

- 1.) Remove all the transportation packaging, except for dollies, from the gantry.

Note: Some sites require floor protection. Locate and install any required floor protection now.

4.2.2 Gantry Prep - For Access Less Than 28"

- 1.) Remove all the transportation packaging, except for dollies, from the gantry.

Note: Some sites require floor protection. Locate and install any required floor protection now.


- 2.) Remove the blue dolly from the left side of the gantry so that the gantry can be positioned closer to the left side wall:

Note: Use Floor Protection for this process.

- a.) Lower the gantry to the floor so that the gantry is resting on the floor.
- b.) Remove the three (3) M14 hex bolts that secure the gantry to the dolly.
- c.) Replace the removed dolly with the shipped black gantry-positioning dolly, and reinstall the three (3) M14 hex bolts.
- d.) Raise the gantry so that it is once again off of the floor.

4.2.3 Gantry Positioning - All Sites


Figure 1-3 Gantry Base Installation Hardware

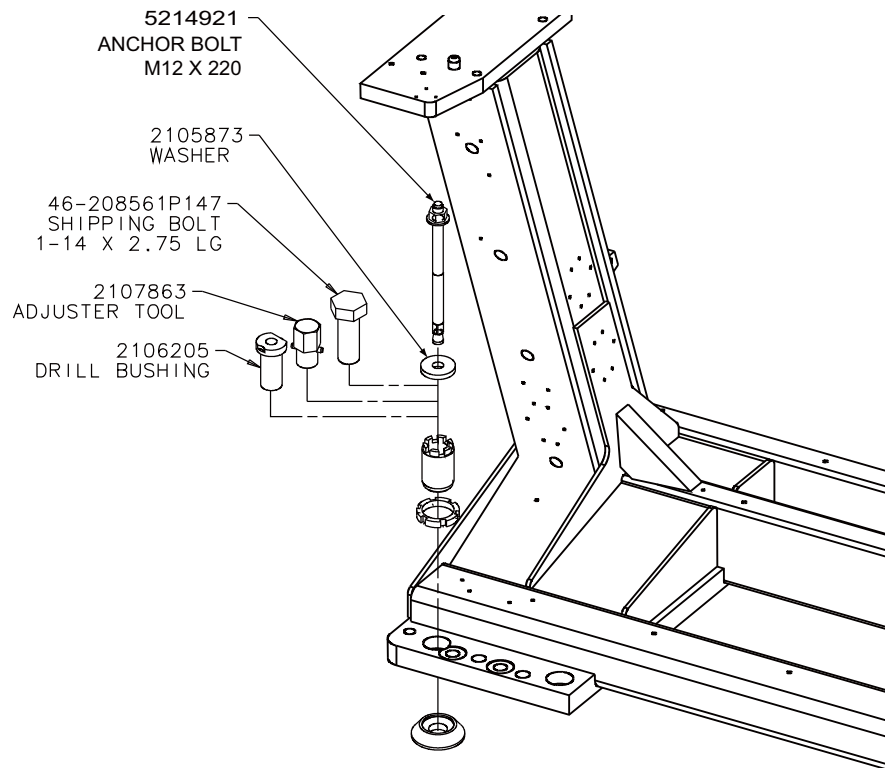

Note: Adjusters are used at each anchor location. Anchor hole ID is 1" (2.5 cm). Void between adjuster and anchor must be filled according to local building codes for seismic application.

- 1.) Position the gantry over the template appropriately.
 - a.) Locate the four (4) leveling pads, and position each of them beneath its associated adjuster.
 - b.) Use the dollies to evenly lower the gantry, until it is just off of the floor (approximately ½" or 12.7 mm). Use a ½" ratchet to raise and lower the dollies.
 - c.) Carefully rotate the gantry into the correct position over the template.

Figure 1-4 Gantry Base Leveling Pads (Starting Positions)

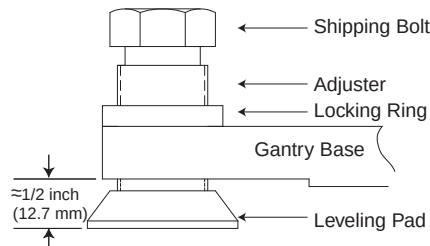


- 2.) Remove the paper templates from the floor and discard properly.
- 3.) Loosen the locking rings and shipping bolts so you can fine-tune the leveling pads to compensate for slight variations in the floor surface.
- 4.) Position the gantry so that the adjusters are centered over their respective holes scribed earlier into the floor.

IMPORTANT: Make certain to route the gantry power cord under the two rear gantry rails, before removing the gantry shipping dollies.

- 5.) Using a ½" ratchet, gently lower the gantry until it rests on the floor, over the marked areas.
- 6.) Using a 14mm hex socket, remove the dollies from the gantry by removing the three dolly bolts found at both ends of the gantry ([Figure 1-5](#)).

Figure 1-5 Gantry Dolly Bolts

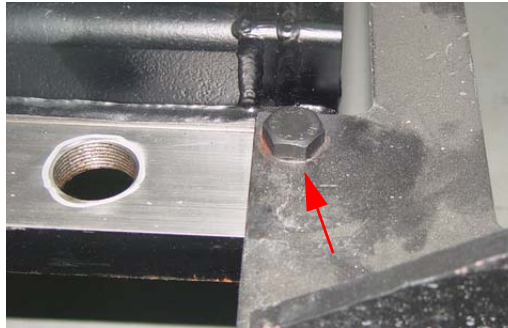


- 7.) Remove the dolly plates on both sides of the gantry. Retain the dolly plates in the service cabinet at the hospital.

Figure 1-6 Dolly Plate



- 8.) Remove the (4) gantry shipping bolts, using a 1½" socket.


Figure 1-7 Gantry Shipping Bolts

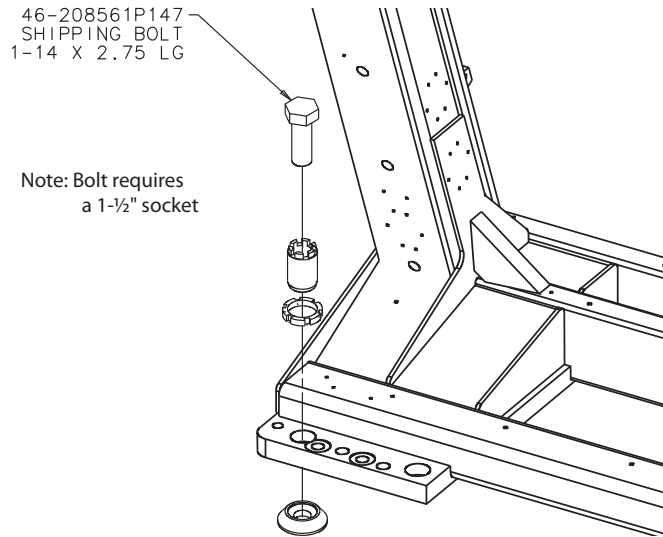

4.3 Level the Gantry

The gantry uses 2 bubble levels that are permanently mounted to machined surfaces on the stationary base to tell when it is level.

Figure 1-8 Gantry Bubble Level

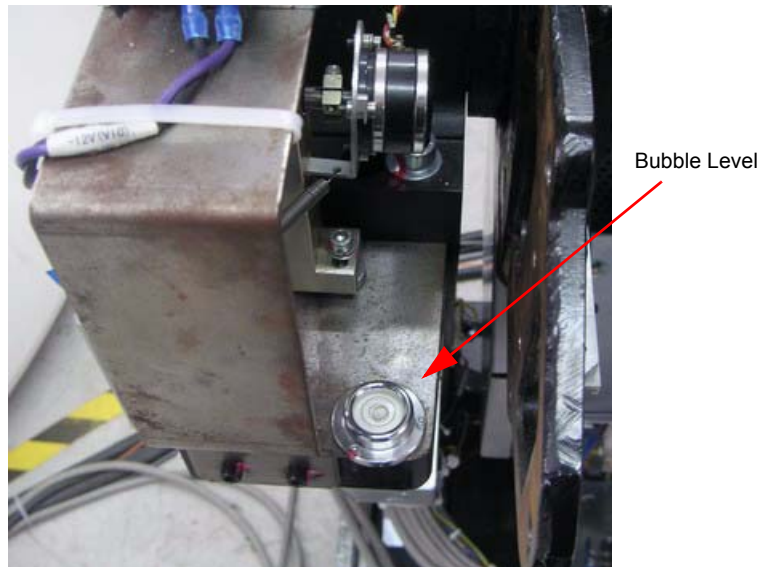


Bubble levels are located on both ends of the gantry stationary base. They're located on the stationary base near a point where the rotating structure pivots mount to the base structure. (See [Figure 1-8.](#)) The gantry is properly leveled when the bubble is centered. (See [Figure 1-10.](#))

- 1.) Loosen all adjuster lock rings (use a spanner wrench or large channel lock pliers).
- 2.) Systematically turn each of the gantry's adjusters (locations 1, 2, 3 and 4 in [Figure 1-9](#)) until both bubble levels are centered left to right, and front to back.
 - Begin by turning each adjuster no more than 1 turn at a time.
 - Use the adjuster tool, 1¹/₈" socket, and the ½" drive ratchet to turn each adjuster. (Refer to [Figure 1-10](#).)

Systematic Procedure for Leveling gantry follows:

- 1.) Level the left side from front to back by turning adjusters #1 and #2.
- 2.) Level the right side from front to back by turning adjusters #3 and #4.
- 3.) Level the side (right or left) that is higher with respect to the other side. Turn both adjusters on a side equally until that side is level. The side should now also be level.


Figure 1-9 Gantry Base “Adjuster” Locations— Top View

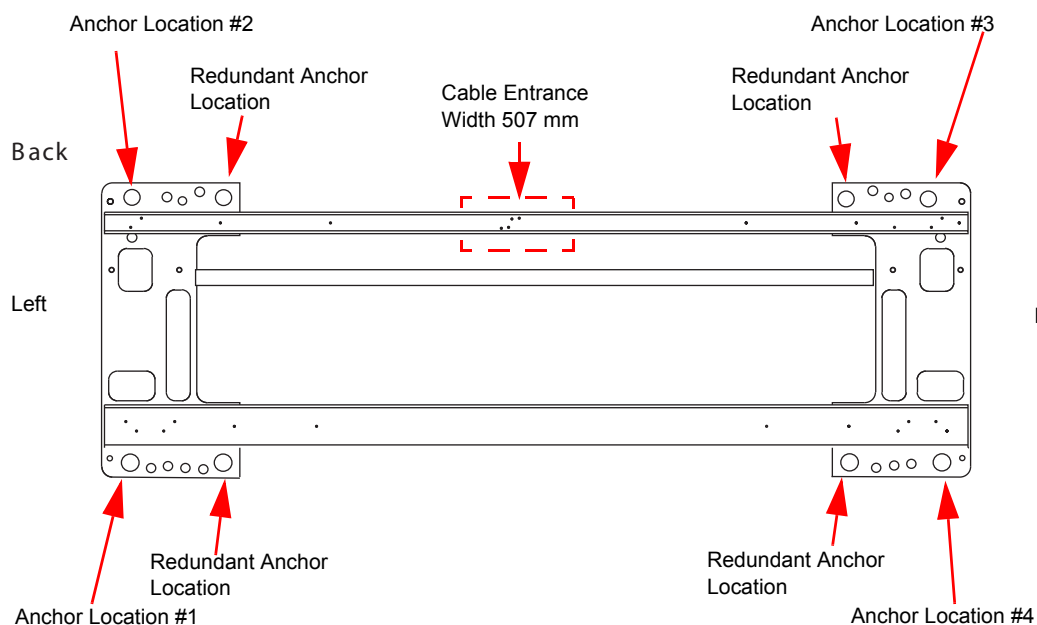

- 3.) When the bubble levels are centered ([Figure 1-10](#)), each of the four (4) leveling pads should be carrying a portion of the gantry weight. Distribution of the gantry weight prevents the base frame from rocking during normal operation. **DO NOT leave any adjuster un-loaded or floating.**

Figure 1-10 Bubble Level Centering

Correct level is 100% of bubble within small circle
Incorrect Level is less than 100% of bubble within small circle

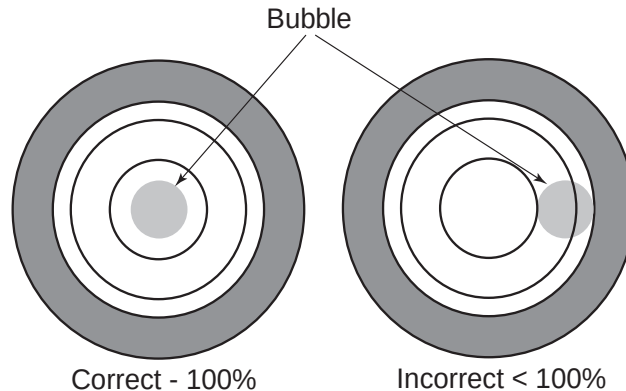

4.4 Gantry Bearing Gap Inspection

All CT systems require a Gantry Bearing Gap inspection before starting electrical calibration. All international gantries are shipped in a wooden shipping crate that should not be removed until it arrives at the installation site. This shipping container is designed to reduce the risk of shipping damage.

4.4.1 Personnel Requirements

REQUIRED PERSONS	PRELIMINARY REQS	PROCEDURE	FINALIZATION
2 (mechanical suppliers or installation team)	15 min	15 min	15 min


4.4.2 Tools and Test Equipment

- Standard tool kit
- Inspection document
- 2.5 mm Allen wrench
- Rear cover dollies (2)
- Flashlight

4.4.3 Damage Indicators

On the inside edge of the black-colored bearing assembly, a mark similar to that shown in [Figure 1-11](#) will be seen, if this is a Rotek bearing.

Figure 1-11 Gantry Bearing - Rotek Label

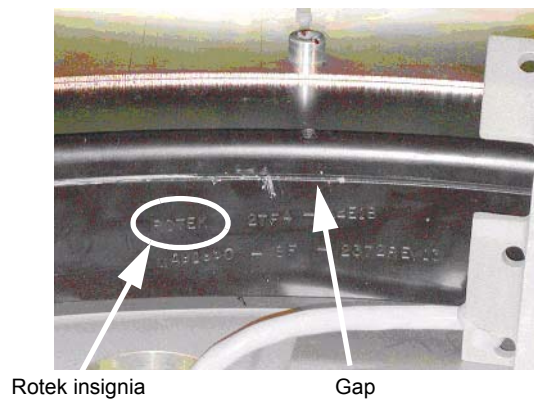


The mark has a serial number in the same format as:

ROTEK 2TF4-44E1B-MA91960-8F-2372-REV13.

The gap to be inspected is shown in [Figure 1-12](#) next to the serial number.

Figure 1-12 Gantry Bearing

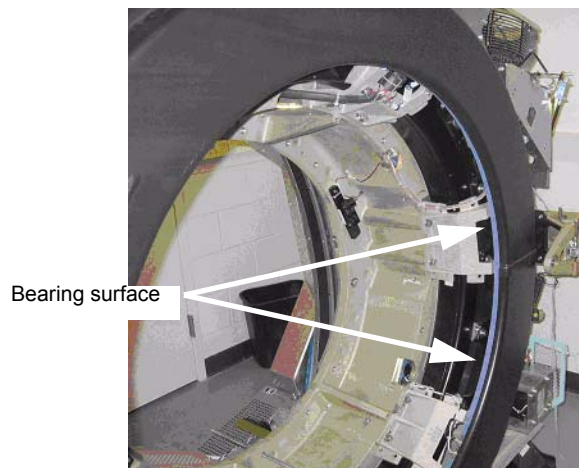


On most systems, a change in the bearing gap does not cause the gantry to make unusual sounds, unless the gap is severe. If the gantry is badly damaged and the gap is severe, it can cause operation issues. Some systems are shipped with shock indicators that must be returned to Milwaukee.

A severe failure may be seen during installation as a problem rotating the gantry.

4.4.4 Porcedure

- 1.) Remove the scan window following the procedure in appendix A, [Section 6.0, on page 128](#)
- 2.) Remove the top and rear gantry covers, following the procedures in Appendix A, [Section 2.0, on page 114](#) and [Section 4.0, on page 124](#)
- 3.) Use a 2.5mm hex wrench as a tool to measure the gap at the positions shown in [Figure 1-13](#). The location of gantry components does not matter. Simply measure four (4) locations 90 degrees apart from each other.

Figure 1-13 Inspection Locations

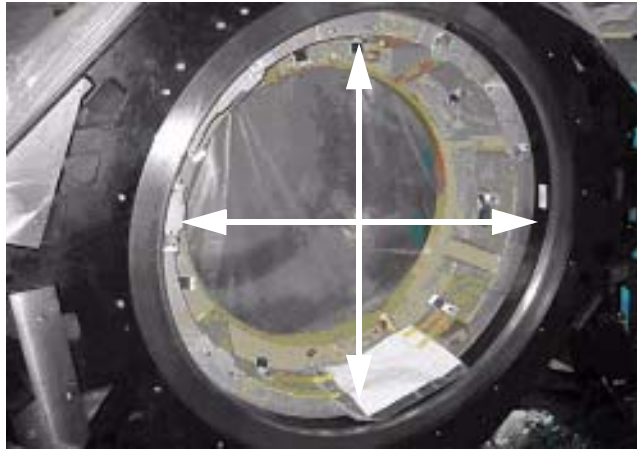

- 4.) If the 2.5mm hex wrench easily fits without effort in the gap, the gap is out of spec. [Figure 1-14](#) shows a gap that is too large. [Figure 1-15](#) shows a gap that is good. Notice that the hex wrench does not fit in the gap in [Figure 1-14](#) but does in [Figure 1-15](#).

Note: Do not use force when putting the wrench in the gap. Either it slips in or it doesn't.

Figure 1-14 Gap Too Large



Figure 1-15 Gap Is Good



- 5.) Replace the top and rear gantry covers, following the procedures in Appendix A, [Section 2.0, on page 114](#) and [Section 4.0, on page 124](#).
- 6.) Replace the scan window, following the procedures in Appendix A, [Section 6.0, on page 128](#).

4.4.5 Finalization - Mechanical Installers

If the Bearing Gap Inspection passes, complete the signoff on the GE Form e4879, Installation Data verification form, that this inspection was completed.

If the Bearing Gap Inspection fails, contact your site FE.

4.4.6 FE Service Action Required

If the Bearing Gap Inspection fails, the mechanical installer notifies the site FE that the inspection failed.

The site FE should:

- 1.) Open a bearing inspection dispatch.
- 2.) Follow the inspection procedure described in this section.
- 3.) Record the bearing inspection results.

If no damage is found, close this dispatch and continue with the electrical calibration procedures.

If the system is damaged, go to the Equipment Delivery Quality web site and follow their instructions.

To enter a damaged in shipping claim, go to this web site:

<http://egems.med.ge.com/edq/home.jsp>

4.4.7 FE Inspection Completion

1.) After the Gantry Bearing Inspection passes, complete the opened service dispatch with the following information:

- Gantry Serial Number
- Gantry Type
- System ID
- Site Name
- Installation date
- Was the Gantry transported to the site in the shipping crate? (Yes/No)
- Was the Gantry lifted or hoisted, were riggers used, or was the Gantry delivered via flatbed wrecker? (Yes/No)
- Number of locations that fail the gap inspection if any: _____


2.) Close the service dispatch.

Should any follow-up be required after this inspection, the site engineer will be contacted directly by CT Engineering.

4.5 Table Prep and Setup


- 1.) Draw a reference line of 564mm position from Gantry Base on the floor as shown in the [Figure 1-16](#).

Figure 1-16 Draw Reference Line

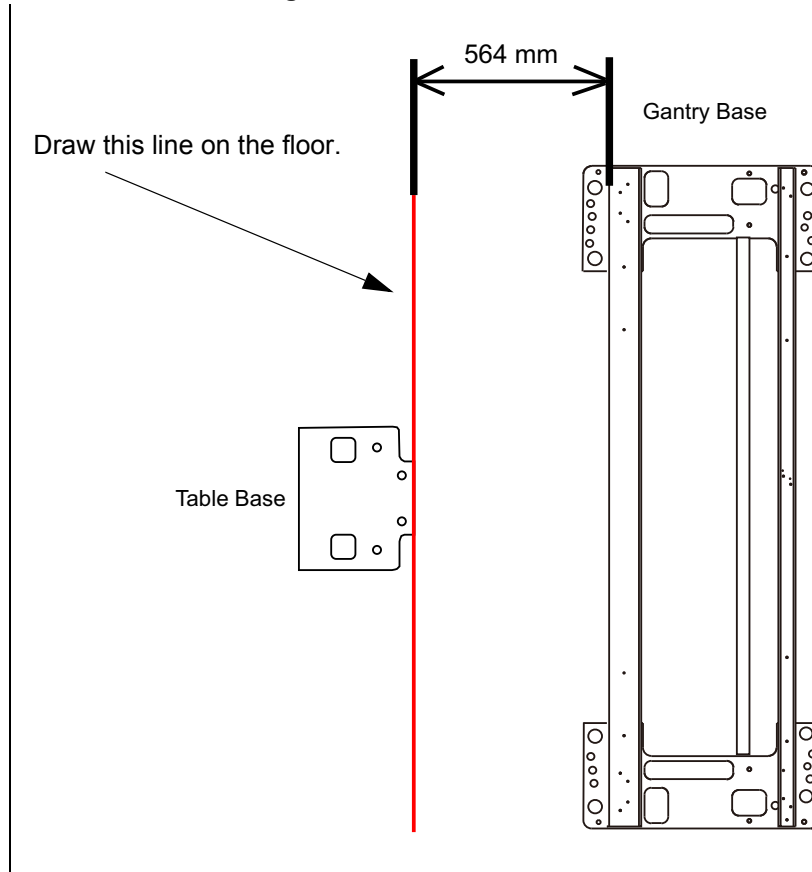


- 2.) Remove all the transportation packaging and boxes, except dollies, from the table.
- 3.) Wheel the table to its approximate position relative to the gantry, using the marks made earlier.
 - a.) Locate the table leveling pads and position them against the base of the table, using the adjusters with a $1\frac{1}{2}$ " socket and $\frac{1}{2}$ " ratchet.
 - b.) Use the dollies to evenly lower the table until it rests on the leveling pads using an $\frac{1}{2}$ " ratchet.

Figure 1-17 Moving the Table

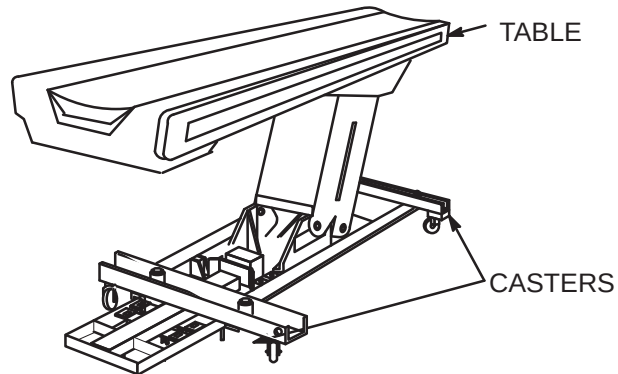


- 4.) Preset leveling pad heights to 20mm, see [Figure 1-18](#).

Figure 1-18 Setting the Adjusters

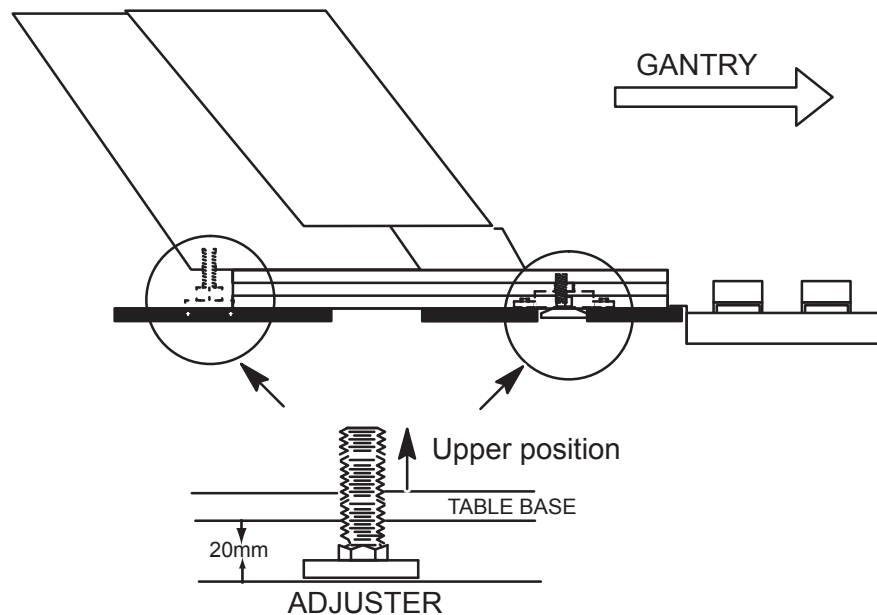

To adjust two leveling pads at rear end of table, please insert the screw driver into the hole located on the bolt then rotate the bolt, refer to the [Figure 1-19](#):

Figure 1-19 Adjust leveling pads

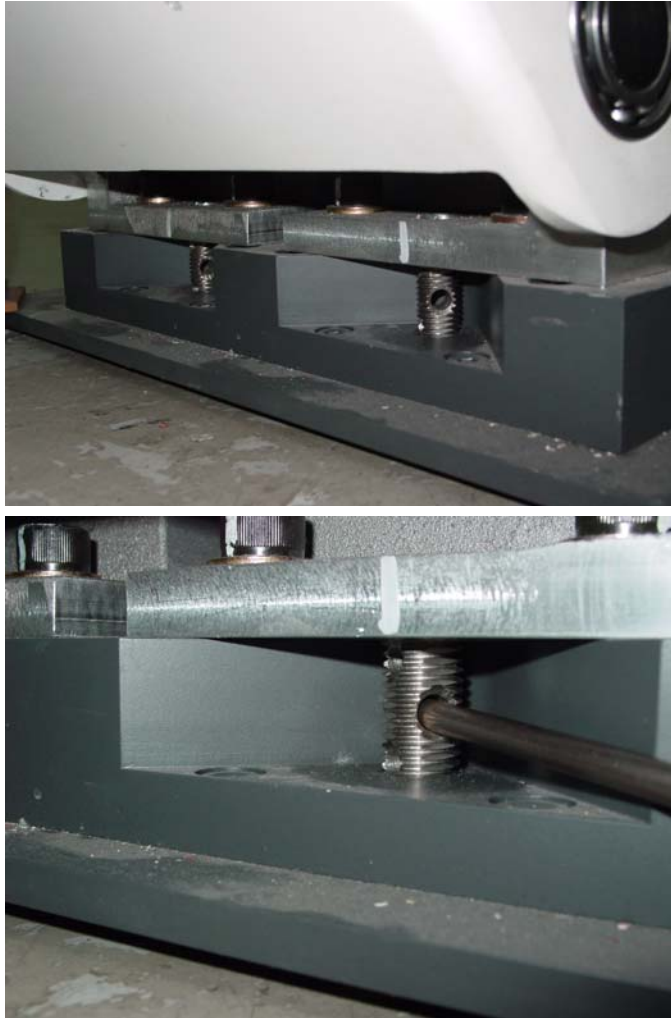

4.6 Cradle Center Procedure


- 1.) Remove a bottom cover under a cradle handle. See [Figure 1-20](#)

Figure 1-20 Remove Bottom Cover



- 2.) Release a latch by pushing a bracket to move the cradle by hands. See [Figure 1-21](#)

Figure 1-21 Release the Cradle

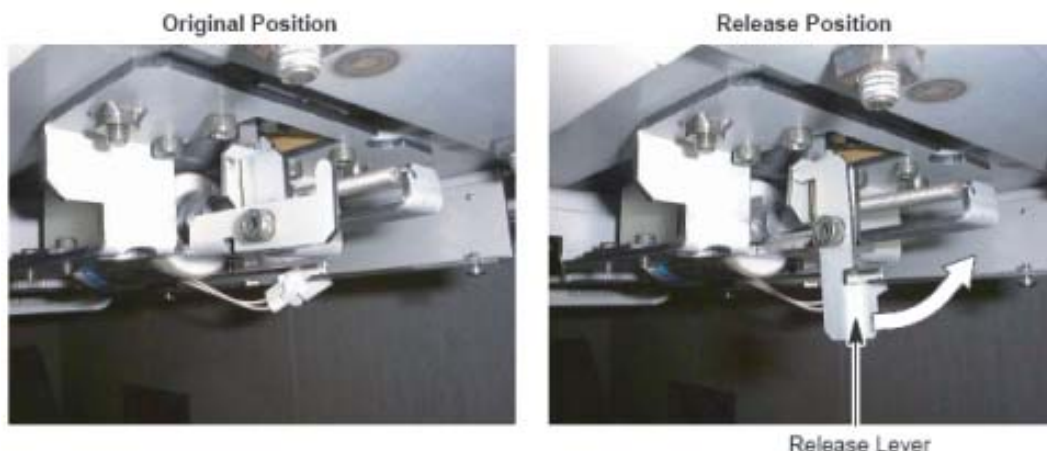

- 3.) Center the cradle on the drive rollers assembly by pushing the cradle into the gantry to its maximum position and back to just before the latch position six times. The cradle should be centered.

You must push and pull a center of cradle hand. See [Figure 1-22](#)

Figure 1-22 Push and Pull Cradle



4.7 Center line on Cradle

- 1.) If your system has no center line mark on the back of cradle, complete the following procedure to make cradle centerline.
- 2.) Using a metric tape measure or other suitable measuring device:
 - a.) Using supplied masking tape, place a 150mm strip on the back of the cradle. See [Figure 1-23](#)


Figure 1-23 Set Strip



- b.) Determine the cradle center-line by measuring the cradle width 428mm (214mm center).

Make two points on the masking tape and make a line. See [Figure 1-24](#)

Figure 1-24 Making the Marker

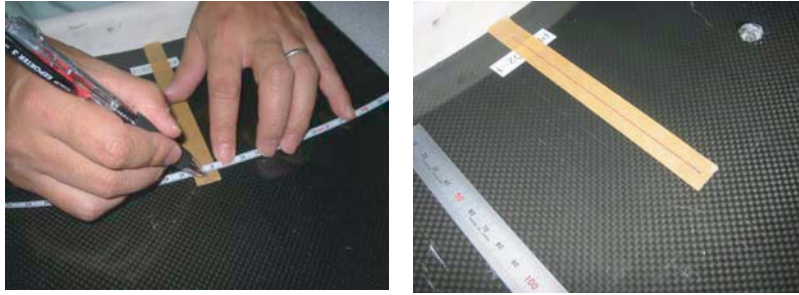

4.8 Install Gantry Alignment Laser and Bracket

NOTICE Use caution while removing the gantry scan window.



- 1.) Rotate the gantry by hand until the collimator face plate is at the 5 o'clock position.

Note: With power OFF, the gantry movement is tight.
DO NOT pin the gantry during this alignment process.

- 2.) With the top and back gantry covers removed, locate the two M10 bolt holes as shown in [Figure 1-25](#). These bolt holes will be used to attach the laser tool to the gantry.
 - a.) The bolts can be installed using a 8mm Allen wrench. Be careful not to bump the alignment light; the mounting space is tight near the alignment light. Tighten until both bolts are snug.
 - b.) Do not drop bolts or the bar on the collimator faceplate. Attach the bar as shown in [Figure 1-25](#).
 - c.) Use a 9" (minimum) level placed on the attached bar to level the bar by rotating the gantry.

Figure 1-25 Alignment Bar Installation Location



CAUTION




Potential for injury.

DO NOT look into the laser.

Use appropriate safety procedures when working with lasers.

- 3.) Attach the laser centering plate onto the laser mounting bar as shown in [Figure 1-26](#). The plate is attached from under the alignment bar using two fixed locators and two thumb screws.

Figure 1-26 Attach Laser Center Plate

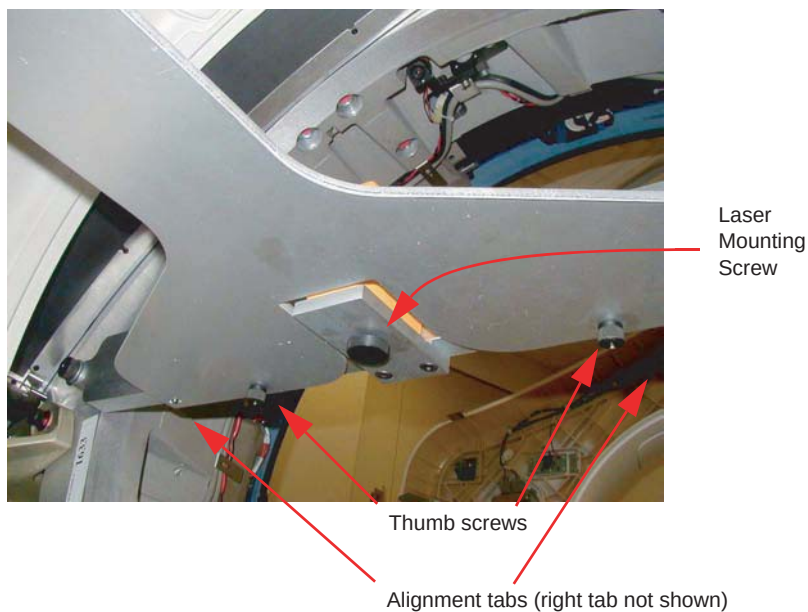


- 4.) When done, insert the laser and turn on the laser using the controls on the back. If the laser is loose when mounted, use a 2 in. piece of Velcro loop (fuzzy) section and attach it to the alignment plate over the attachment screw. Remount the laser and it should fit snugly without moving.

- 5.) When pressed, the **ON** button steps through four different beam profiles and “Self-Leveling Off”. Press the **ON** button until the “|” beam shows. It will be used for this operation.


Times pressed	Function	Notes
1	—	Self-leveling on
2		
3	+	
4	Self-leveling off	Do not use

Figure 1-27 Laser On Button

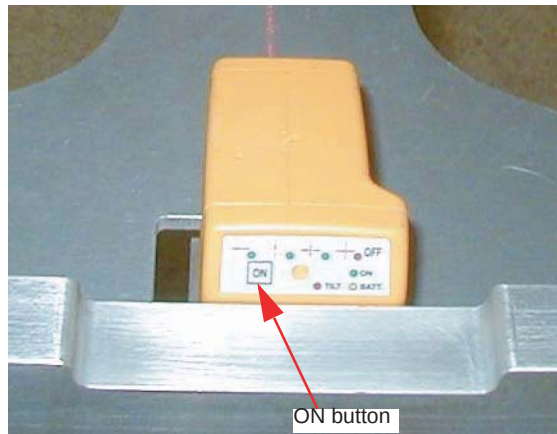

- 6.) Align the laser by carefully rotating the laser base assembly so that the “I” beam shines through the alignment sight mounted on the end of the alignment plate. Use the locking screw on the bottom of the alignment bar to secure the laser to the bar, as shown in [Figure 1-28](#).

Figure 1-28 Laser Centering

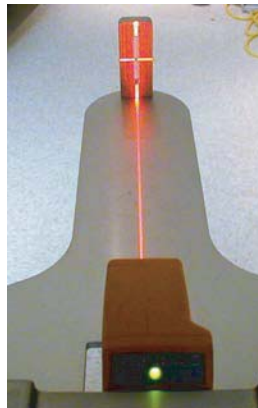


Note: When tightening, the laser may move. Use caution to prevent any movement, as this can result in drilling the table anchor holes in the wrong location.

- 7.) After the laser is centered, notice that the laser beam also appears on the back wall. Place a piece of masking on the wall and carefully mark a line on the tape where the laser appears. This line is later used in the table alignment. This line is also useful in determining if the laser unit moves during the alignment process.
- 8.) Remove the alignment centering plate and store it in the alignment case.

- 9.) Using a chalk line, mark a table center line on the floor along the laser light shining on the floor.
- 10.) Check the reference line for table 564mm Z-distance, as shown in [Figure 1-29](#)

Figure 1-29 Checking the reference line

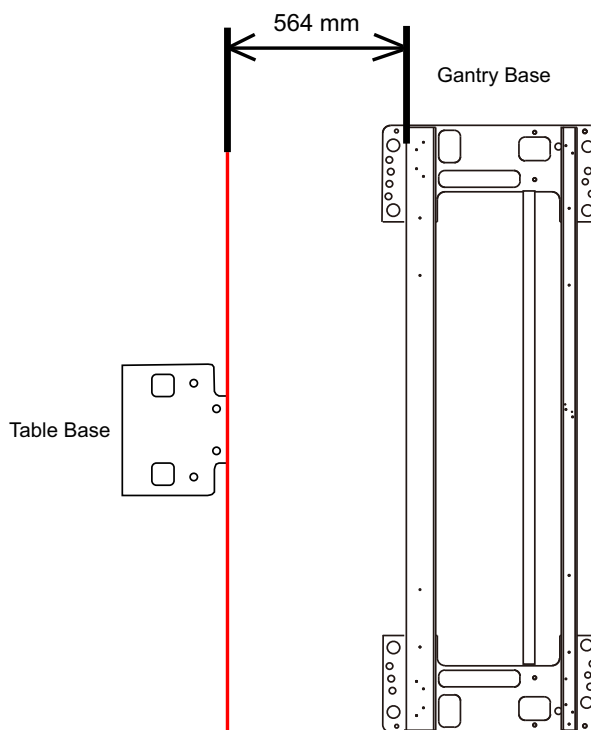



- 11.) Turn off the laser, but do not remove.

4.9 Install the Table Cradle Laser Alignment Plates

Required tool: TG-alignment tool 5180876

- 1.) Install the front table laser alignment which has a cross slot to the front of the cradle.
 - a.) Loose adjustment knob to make the slot plate free.
 - b.) The slot in the front of the tool should be over the cradle center-line. See [Figure 1-30](#)

Figure 1-30 Front Table Laser Alignment

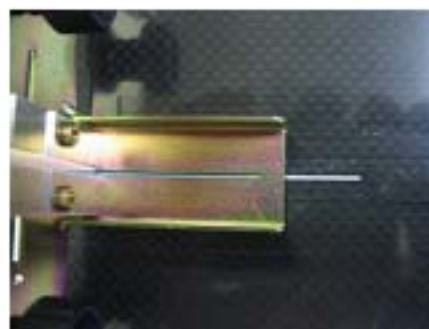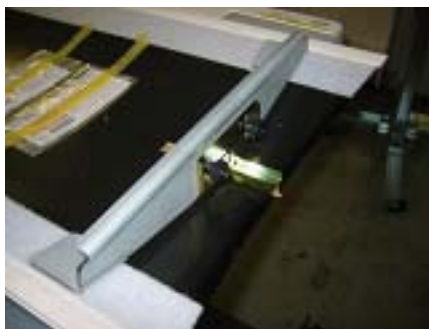

- c.) Tighten adjustment knob.

- 2.) Install the rear table laser alignment which has a cross target to the back of the cradle.
 - a.) Loose adjustment knob to make the slot plate free.
 - b.) The slot of the tool should be over the cradle center-line. See [Figure 1-31](#)

Figure 1-31 Rear Table Laser Alignment



- c.) Tighten adjustment knob.
- 3.) Check that table base is centered over the table center line, and the base is on the 564mm line made on the floor, see [Figure 1-29](#).

4.10 Level and Center the Table to the Gantry

4.10.1 Required Tools

- 2' Carpenter's Square
- 2' Level
- Masking Tape
- Flat-blade Screwdriver
- Craftsman Laser Trac™ Level with Carrying Case and Laser Enhancing Glasses
Sears item #00948247000
Mfr. model #48247
\$39.99
- 4' Level
- 9" Level
- Metric Tape Measure
- 1' Pry Bar

4.10.2 Alignment Conditions

- Before you start, turn on the laser and check that the beam is still on the mark placed on the wall. If not, reset the laser.
- If the mark is not present, use a measuring tape and place a 100 mm piece of masking tape on the cradle at the 1000 mm and on the laser line.
- The gantry must be at zero degrees.

4.10.3 Alignment Specification

Mechanical alignment must be within ± 3 mm
Table cradle travel must be parallel
Table must be perpendicular to the gantry

4.10.4 Procedures

NOTICE Avoid leaning on the cradle during this procedure.
DO NOT pin the gantry during this alignment process.
This procedure is for systems mounted on 102 mm (4 in.) concrete floors only!

Note: If the floor covering was not properly removed with the glue removed, or the levelers were not centered over the floor cutouts, the leveler may become trapped against the edge of the floor covering, causing the table to become unlevel. If this happens, move the table and enlarge the 102 mm (4 in.) floor cutout for the table. Glue removal is important and aids in moving the table to its final location in accordance with the floor levelness specification.


- 1.) Place the bubble level on the table base . See [Figure 1-32](#)

Figure 1-32 Place Bubble Level



- 2.) Turn on the laser's "I" beam (vertical beam) by pressing the ON button 2 times
- 3.) With the laser on, the laser light should shine through the front tool down the cradle connecting the front and rear center lines and onto the center of the rear alignment tool.
- 4.) Raise or lower the table as needed using the two front and two rear adjusters so that this line is shining 100% in the machined grooves of each alignment tool.(it is easier to raise, move, and lower the table using the table dolly. This prevents the adjuster levelers from gripping on the floor adhesive, making it difficult to move.)
 - 1.) The light will shine vertically down the table.
 - 2.) Both center lines on the cradle front and rear should be connected when done.
 - 3.) Level of the table influences a center position, The vertical laser alignment light must match and be centered on the white line on the front of the table and in the groove of back alignment plate. See [Figure 1-33](#)

Figure 1-33 Front and Rear Adjuster

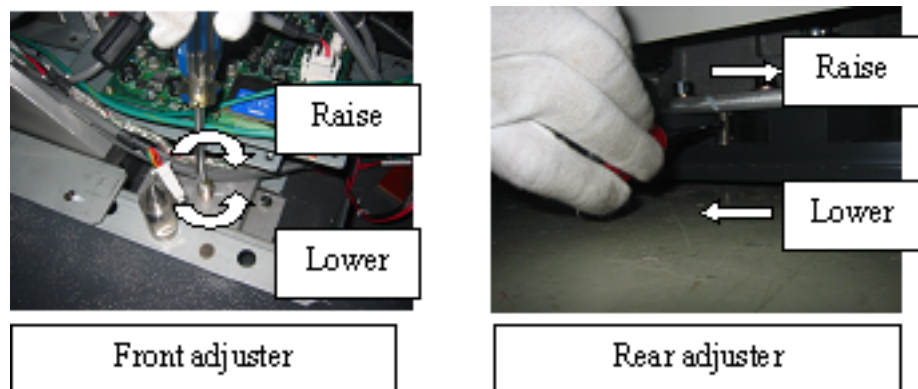

- 5.) Be sure that all adjusters bear part of table's load.

4.11 Drill the Anchor Holes

WARNING POTENTIAL FOR PATIENT INJURY.



**IMPROPERLY SECURED TABLE MAY TIP, DISLODGING PATIENT.
PROPER ANCHORING IS KEY TO MAINTAINING PATIENT SAFETY DURING
SYSTEM OPERATION.**

4.11.1 Notes to Mechanical Installers

GE-provided floor anchors are designed for use **ONLY** on concrete floors that meet the 4-inch concrete floor requirement. Supplied floor anchors must be installed by a trained contractor, and shall be set to a minimum depth of 3 inches at each anchor point. ANY anchors having more than 1-inch of thread showing above the nut, when torque is set to 55 lb.-ft, shall have a second anchor installed in the closest adjacent hole. This is because the minimum anchor engagement length in the concrete was not met. The second anchor shall be installed to the standard depth and torque specification. **Do not cut anchor bolts that extend longer than the 1-inch limit.**

4.11.1.1 Note 1: Alternate Anchoring

If at least four anchors cannot be set for the gantry, and at least four anchors for the table, then the installer must inform GE that the minimum anchoring cannot be met. Additionally, a structural engineering contractor must be engaged to determine the anchoring method and to certify that their anchoring meets the stated GE minimum load requirement and torque specification.

4.11.1.2 Note 2: Non-Concrete Floors

All other anchoring methods — on floor types other than the concrete minimum — must be determined at the customer's expense, by a structural engineering contractor. The anchoring method must be certified to meet the stated GE minimum load requirement and torque specification.


4.11.1.3 Note 3: GE Notification

It is not the role of the mechanical contractor or installer to determine acceptable methods to install or anchor equipment on non-4-inch concrete floors. The appropriate GE contact person shall be notified that the facility's floor type **DOES NOT MEET** the installation mounting requirement for the installation procedure (described in this Installation Manual), and therefore the table-gantry mounting process **CANNOT** continue.

4.11.2 Gantry Drilling

- 1.) Make certain that the gantry scan window is in place, to prevent dust from entering the gantry.

Figure 1-34 Gantry Scan Window

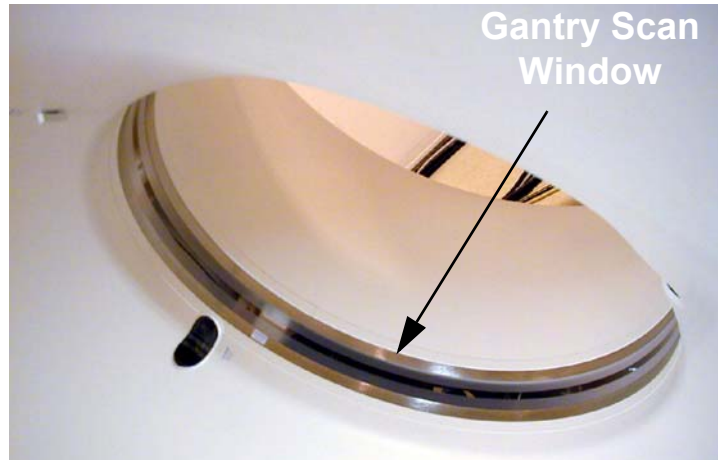


- 2.) Use a piece of tape to mark the drill bit depth of 220 mm from the tip of the 12mm masonry drill bit.

NOTICE
Potential for
Equipment
Damage from
Dust

You must cover all electronic assemblies in the table base prior to drilling to prevent damage due to the dust created during drilling.

- 3.) Use the 12mm drill bit to drill all four (4) anchor holes (refer to [Figure 1-36](#) for locations.)
 - Review [Figure 1-35](#), prior to drilling.

Figure 1-35 Drilling Position

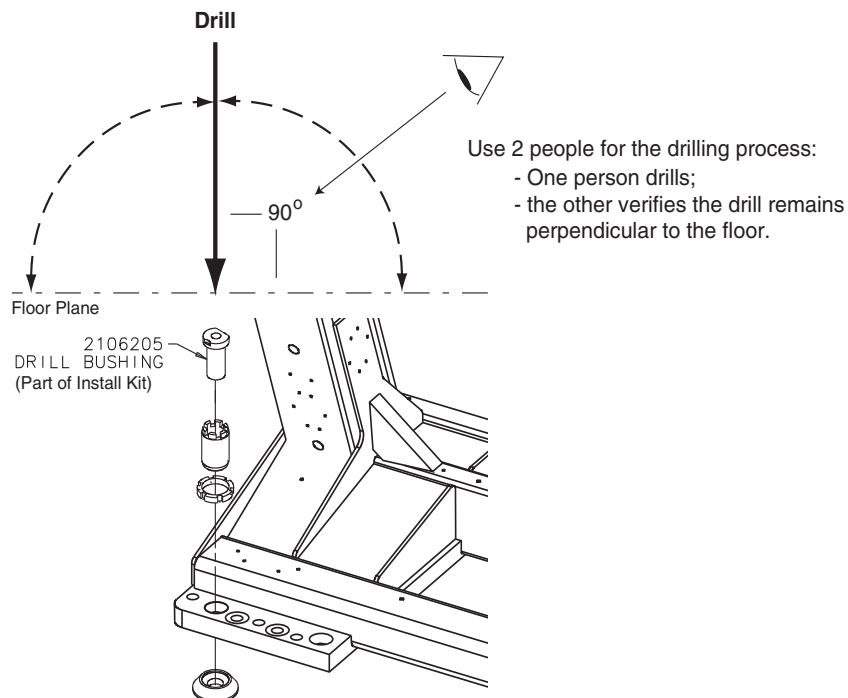

- **Important** - Drill the holes perpendicular to the floor. While one person drills the holes, position a second person to watch the relationship between the drill bit and floor. Make

sure the bit remains absolutely perpendicular to the floor throughout the drilling operation. Always use the mechanical guide supplied with the alignment kit.


- Stop drilling every 15 or 20 seconds and clear the hole of debris. This helps to prevent binding of the drill bit.
- a.) Drill each hole until the mark on the drill bit is even with the top of the drill bushing. (195mm)
- b.) Place the drill bushing inside each adjuster, to keep the hole vertical, and centered within the adjuster.
 - * Take care not to injure yourself on the gantry cover brackets.
- c.) Stop the drill every 15 or 20 seconds and clear the hole of debris to prevent binding of the drill bit.

- 4.) Vacuum all debris from the inside of the hole:

Note: A drywall dust filter should be used on the vacuum.

- a.) Place the funnel tip inside the hole; place the vacuum hose in the funnel.
- b.) Continue to vacuum while you drill, if you can, to keep gantry and table as free of dust contamination as possible.
- c.) When you finish clearing the anchor hole, vacuum the debris from the surrounding area.
- d.) Stop drilling frequently to let the drill bit cool.
- e.) All holes must be a minimum of 100 mm deep (see [Figure 1-41, on page 56](#)). Recheck all holes.

Figure 1-36 Anchor Locations, Standard

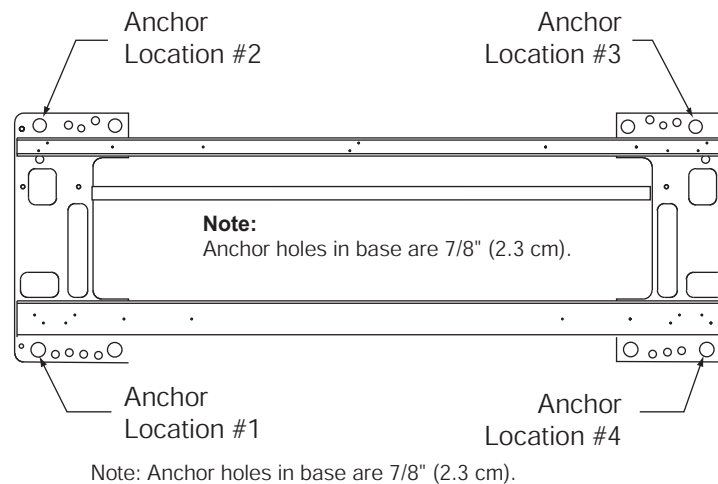

4.11.3 Table Drilling and Anchoring

When the room-layout template is used, the anchor bolts have been installed. If it is, perform the following. If it is NOT, perform all procedures of this section.

- Table Leveling (This section)
- Tighten the anchor nuts with the insulating plate and washer.

TORQUE: Approx. 5 N-m

NOTICE

Do not tighten the anchor nuts to the specified torque (50 N-m). To over-tighten nuts will damage the frame of the Gantry and also bring an out-of-level Gantry.

- Verify that the Table to Gantry Alignment should be correct.


This anchoring procedure is one for the anchor kit shipped with the system.

NOTICE The minimum Tensile load strength of the anchor must be 13000 N.



- 1.) Verify that the table flapper is removed to make space for anchoring. If it is not, remove the table flapper by removing the four flapper Assy mounting bolts.

Figure 1-37

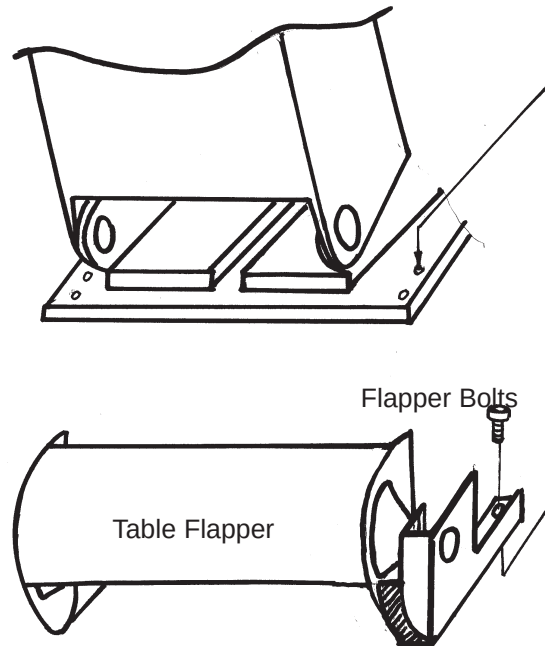


- 2.) Prepare and tape 12 mm concrete drill bit (2203081: HILTI, TE-CX 12/17 or equivalent).(Refer to step 5.), on page 53)
- 3.) Make holes at the anchor position as shown in Figure 1-40, using electric drill with the bit prepared in the previous step, so that the lower edge of the tape just touch the floor surface.

NOTICE During drilling, always keep the drill bit perpendicular to the floor. With one person drilling the hole, a second person can visually insure perpendicularity throughout the drilling operation. Clear the hole of debris several times while drilling to prevent the drill bit from binding.




- 4.) Attach a piece of tape to the 12 mm concrete drill bit (2203081: HILTI, TE@CCX 12/17 or equivalent) in order to assure a minimum of 100 mm depth of anchor hole.
- 5.) 12 mm Concrete Drill Bit Preparation

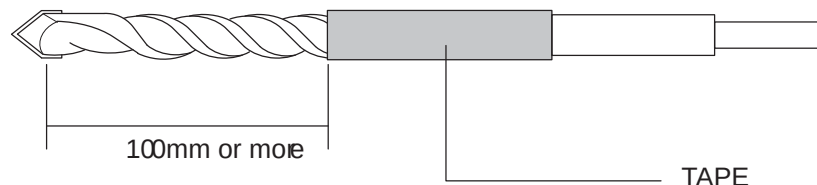

- 6.) Use the vacuum cleaner to suck up the concrete dust from the hole. Also vacuum dust from the surrounding area when finished.
- 7.) Place the bubble level at the positions.
- 8.) Level the Table using the four adjusters. (The front sides of the adjuster contain the lock nut. After level adjustment, tighten the lock nuts securely.)

Level Specification: ± 0.5 mm/1000 mm

Figure 1-38 Table Adjuster and Bubble Level Position

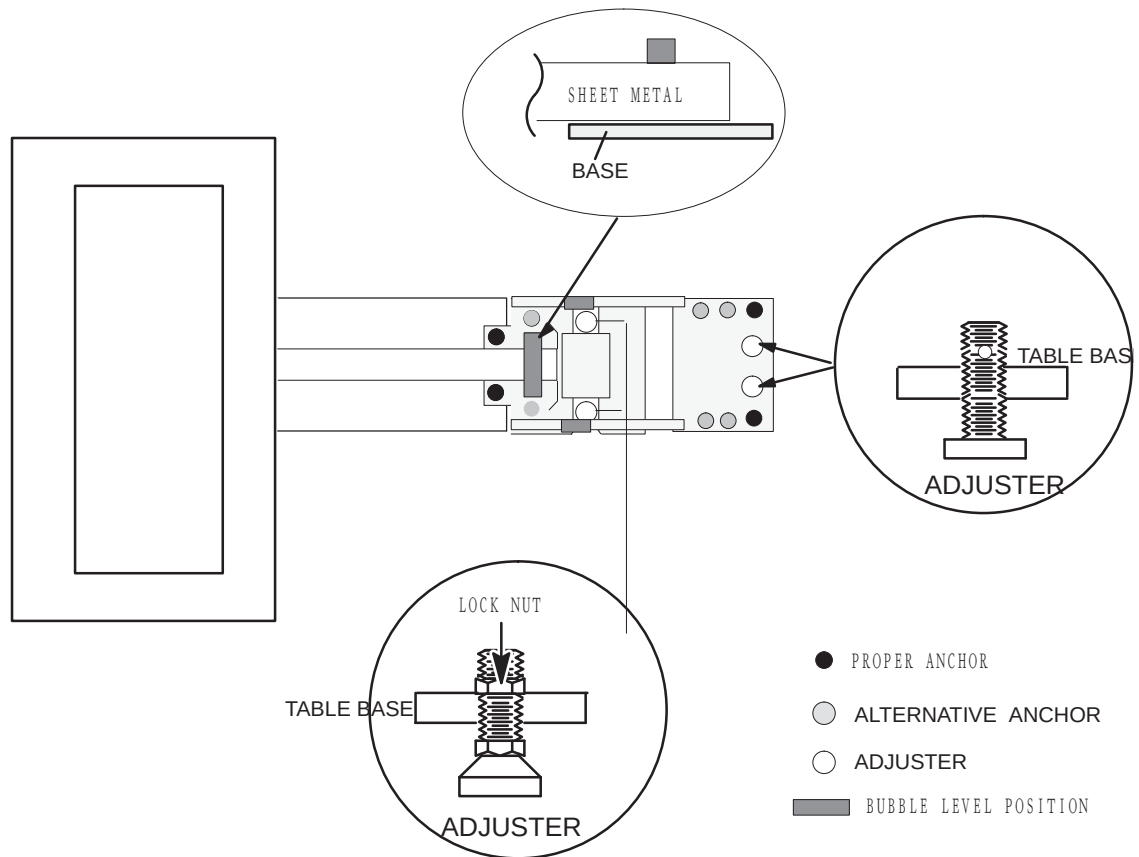

9.) Check that the Cradle is also leveled.


4.12 Gantry & Table Redundant Anchor Holes

If you cannot use one of the adjuster anchor holes due to structural interference, such as reinforcement bars in the concrete, you must use one of the redundant anchor locations, as shown in [Figure 1-39](#).

- The gantry requires a minimum of four (4) anchors, one (1) in each corner of gantry base.
- The table requires a minimum of four (4) anchors, one (1) in each corner of table base.

If you must use a redundant anchor hole in the gantry, you must remove the gantry covers. See [Appendix A](#), for gantry cover removal.


Figure 1-39 Table Base Redundant Anchor Locations

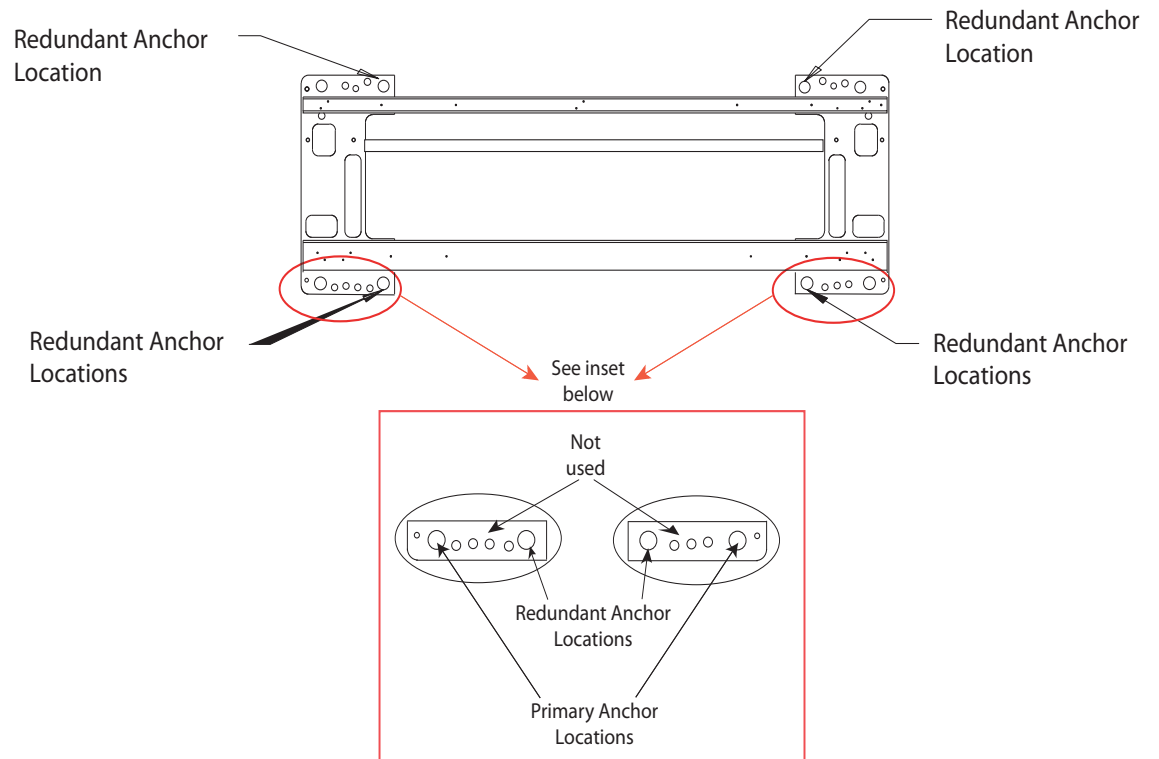

Mark the 12mm drill bit with a piece of tape 170 mm ($6\frac{11}{16}$ ") from the tip of the bit, and follow the instructions in [Section 4.11, on page 50](#).

Note: If any obstruction is found under the proper anchoring position, stop drilling and use the alternate position (see [Figure 1-40](#)). Use the alternate positions only when the case its corresponding proper position is not possible.


Figure 1-40 Table Anchor Positions (Proper and Alternative)

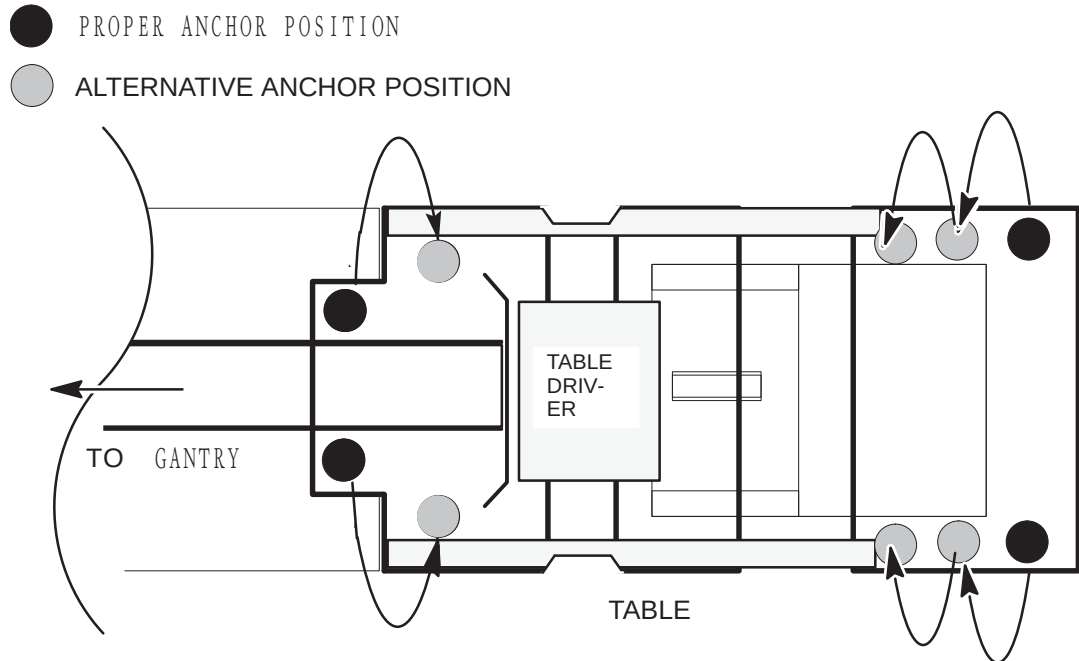

4.13 Gantry Anchoring

Recommended - Use "HSA - KA" anchors P/N 5214921 (12mm dia. by 220mm long) as shipped with the BrightSpeed Select for this procedure.


Figure 1-41 Gantry Base Anchor Assembly

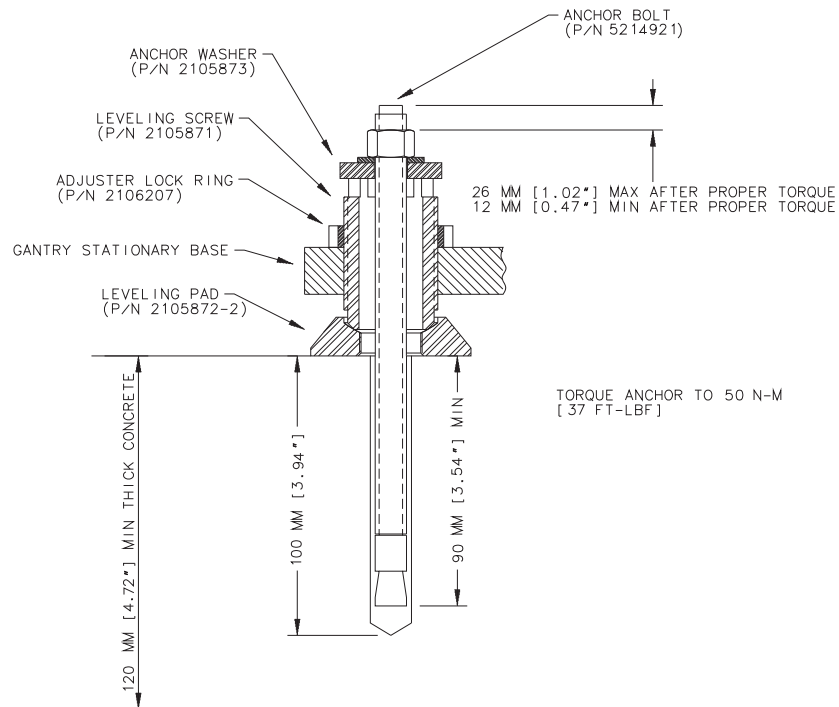

- 1.) Remove the alignment tool and repack the install kit.

- 2.) Assemble the anchors before you install them. Refer to [Figure 1-41](#).
 - a.) Remove the nut and washer from the anchor.
 - b.) Add a 1/4" thick washer (PN 2105873) under the regular anchor washer.
 - c.) Reassemble the anchor washer and nut and position nut so top is flush with threads of anchor before hammering into hole.
- 3.) Fixed site Installations must be centered under mounting holes in table base, with 3mm of space visible surrounding the tapped bolt hole. [Figure 1-42](#).

Figure 1-42 Center tapped holes under mounting holes in table base

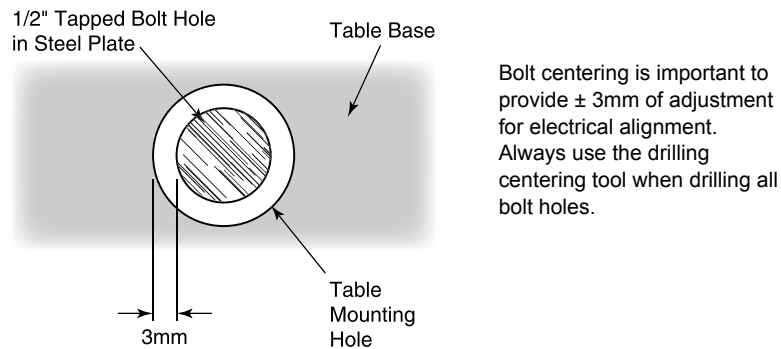

4.14 Alignment Recheck

Note: Alignment is critical. Recheck carefully.

- 1.) Turn on the alignment tool and recheck alignments. The table alignment must be the same as in [Procedures, on page 49](#). If re-leveling is required, repeat this procedure. Using the bubble levels, make adjustments as required to maintain required alignment.
- 2.) Once alignment has been verified, torque all mounting bolts. Tighten anchors and torque to $50 \pm 6 \text{ N-m}$ ($37 \pm 5 \text{ ft.-lb.}$)
- 3.) Remove the laser tools.
- 4.) Reinstall all the removed table panels and hardware.
- 5.) Reinstall the gantry rear cover.



Note: If you cannot replace the lower table cover because the floor interferes, adjust all of the table and gantry levelers by half-turn increments to raise the table/gantry until the lower table covers clear the floor. Then return to the alignment sections to level the gantry, level the table, and tighten the locking rings, respectively.

Section 5.0

Rear Entry Cable Box

A rear entry cable box (B7850RC) is used when the cables to the gantry cannot be brought-up inside the gantry base. The box is not supplied with the system and must be ordered separately. The rear entry cable box frame attaches to the gantry base. See [Figure 1-43](#).

Figure 1-43 Rear Entry Cable Box

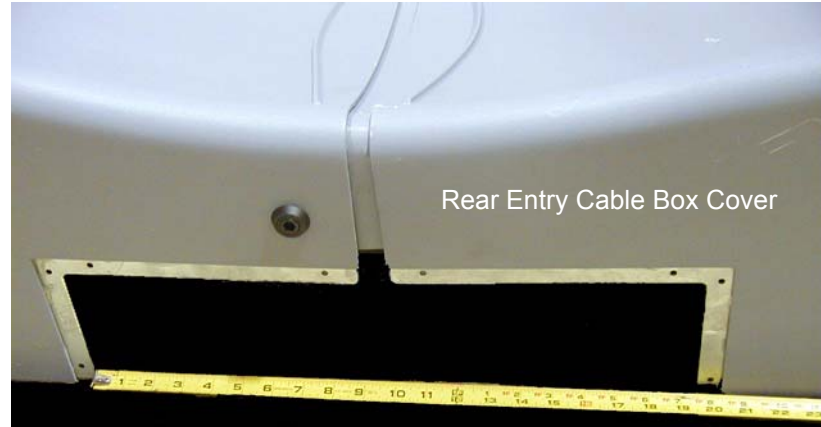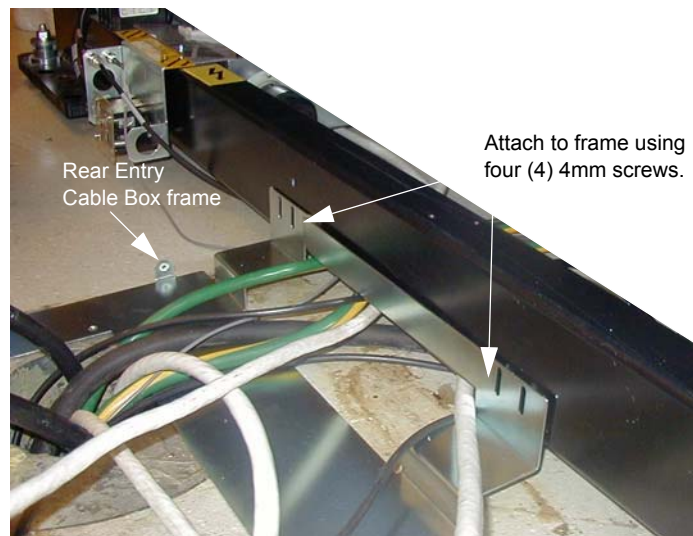

The filler plate may need to be cut in order to fit the rear cable box into the cover opening. Use the supplied white molding around the filler plate.

Section 6.0

Install Table Footswitch Assembly

Install the table foot-switch assembly as shown in [Figure 1-44](#), along with filler strip.

- 1.) Connect foot-switch cables.

Figure 1-44 Install Table Foot-switch



- 2.) Install front bottom cover and footswitch pedals.

Figure 1-45 Install Front Bottom Cover




Figure 1-46 Install Footswitch Pedals



Section 7.0

Remove Gantry Tilt Bracket

Figure 1-47 Gantry Tilt Bracket Removal

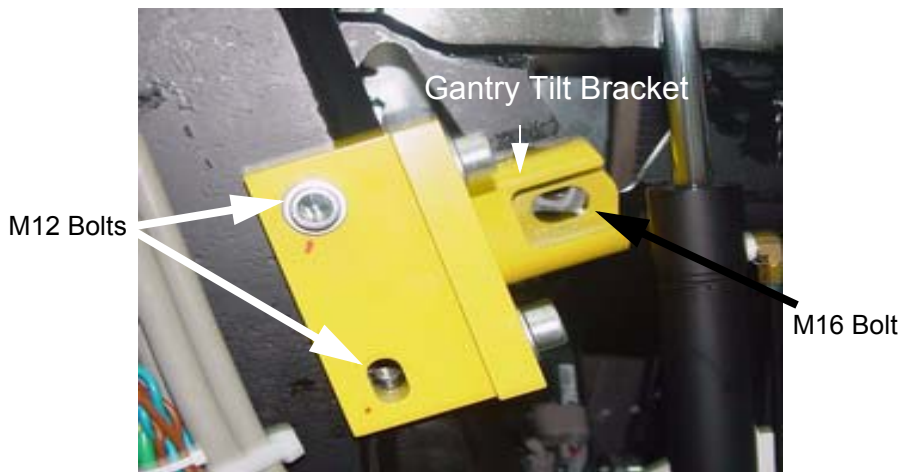

- 1.) Refer to [Figure 1-47](#). Remove the M12 bolts using a 12 mm Hex wrench.
- 2.) Loosen the M16 bolt 1-2 turns and check the Gantry tilt bracket, it should be loose to the touch. If loose continue with step 4.

CAUTION



If tilt bracket is not loose, stop and put the M12 bolts back in and tighten tilt bracket back in place. If there is a load on the tilt bracket, removal may cause the gantry to suddenly tilt all the way back due to a possible lack of hydraulic pressure.

- 3.) Check the hydraulic connections for leaks or lack of fluid. You will have to wait until the system can be energized to use the tilt controls to relieve the load on the tilt bracket prior to removal. Do not use force to remove the bracket.
- 4.) If the bracket feels loose, remove the M16 bolt using a 16 mm Hex wrench.
- 5.) Remove the bracket.
- 6.) Close the gantry covers and reinstall the scan window.
- 7.) Store brackets in the gantry base.

Section 8.0

Position the Power Distribution Unit

WARNING



LOCKOUT/TAGOUT IS REQUIRED BEFORE PERFORMING THIS TASK. USE THE SUPPLIED LOTO KIT.

- 1.) Roll the PDU into position on its permanently mounted casters. Leave at least 15.5 cm (6") between the PDU and back wall to allow cooling air to circulate.

Note: Connecting the primary incoming power (steps 2 through 5, below) is performed by the customer's electrical contractor.

Table 1-3 Contractor Connections

Connection or Wall Box	AWG #	Connection From	Connection To PDU	Installed & Checked
TS1	#1	PDB-A	TS1-1	
	#1	PDB-B	TS1-2	
	#1	PDB-C	TS1-3	
	#1/0	GND	N/G (Do NOT connect anything to neutral point.)	


- 2.) Run the main input power conductors and ground through flexible metal conduit (attached between the PDU chassis and room duct-work) so you can move the PDU away from the wall during service.

Figure 1-48 Flexible Conduit for PDU Power



- 3.) Locate the hole cover plate in Box 1 and attach the flexible metal conduit to the PDU.
- 4.) Cut the three phase and 1/0 ground wire to size.
- 5.) Attach cables as shown in [Figure 1-48](#) and [Figure 1-49](#).

Figure 1-49 PDU Area Locations

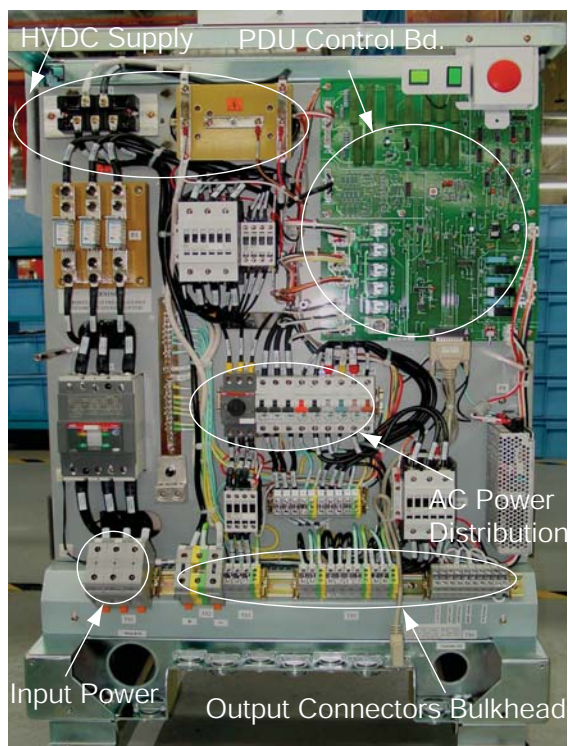

Section 9.0

Install Operator Console Without FWS Table

9.1 Unpack Console

- 1.) Remove all items from the console.
- 2.) Remove all packing materials and discard.

Note: For the Operator Console shipped without skid package, go to [9.2](#).

Figure 1-50 Console boxed on skid



- 3.) Place the step-board under the front edge of the skid and step on it to raise the front edge of the skid as in [Figure 1-51](#).

Figure 1-51 Step-board used to raise front edge of skid



- 4.) Remove the two front cushions from the bottom of the skid. Refer to [Figure 1-52](#)

Figure 1-52 Cushion on bottom of skid



- 5.) Seismic Brackets are shipped with system. Remove the Seismic Brackets from each side of the console. Refer to [Figure 1-53](#). Save brackets for use later if you need to mount the console to the floor.

Figure 1-53 Removing the Seismic Brackets




- 6.) Lift up on the strap on the front of the step-board ([Figure 1-51](#)) to lower the skid. Remove the step-board.
- 7.) Ensure the console stabilizers are in line with the notched portion at the front of the skid. This will allow enough clearance to smoothly roll the console down the ramps.
- 8.) Move console to installation location.
- 9.) Install the anti-tip legs of the console: Mount the two anti-tip legs on the console: 2318681-2 (left), 2318681-3 (right), secure three hex screws from console side.

Figure 1-54 Mount anti-tip leg on the console



- 10.) Attach the two pads to the console stabilizers. Refer to [Figure 1-55](#).


Figure 1-55 Console Stabilizer Pads

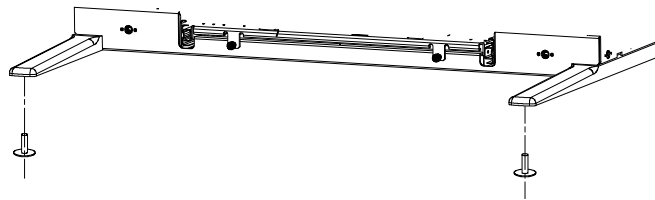

- 11.) Screw the pads all the way into the legs

9.2 Remove Console Covers

Figure 1-56 Front Cover Removal - Global Consoles

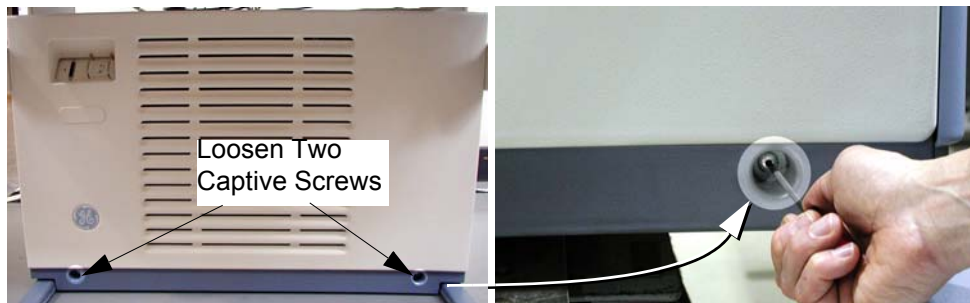

- 1.) Loosen two captive screws at bottom of console.
- 2.) Rotate bottom of cover outward and upward until it may be lifted free of the console at the top.

9.3 Adjust Table Top Height and Position Console

- 1.) Refer to [Figure 1-57](#). The console normally arrives with the bottom hole on the monitor top aligned with the fourth bolt from the bottom.


Figure 1-57 Global Console tabletop adjustments



- 2.) For optimum operator comfort, align the bottom hole of the keyboard table top with the fifth hole from the bottom on the console.
Your site may have different requirements and you may have to adjust the monitor top and/or keyboard table top up or down from this position.
 - Always select a table top height that permits the operator to see the patient on the table.
 - Keep the one bolt-hole relationship between the monitor top and the keyboard table top. Fasten the keyboard table top one (1) hole lower than the console table top (Figure 1-57).
- 3.) Install the console side covers.
Side covers are left and right specific.
- 4.) Move the console into its final position, maintaining a 6" (15 cm) minimum distance between the console and the wall.

9.4 Attach Keyboard Table Top

Figure 1-58 Keyboard Table Top Placement

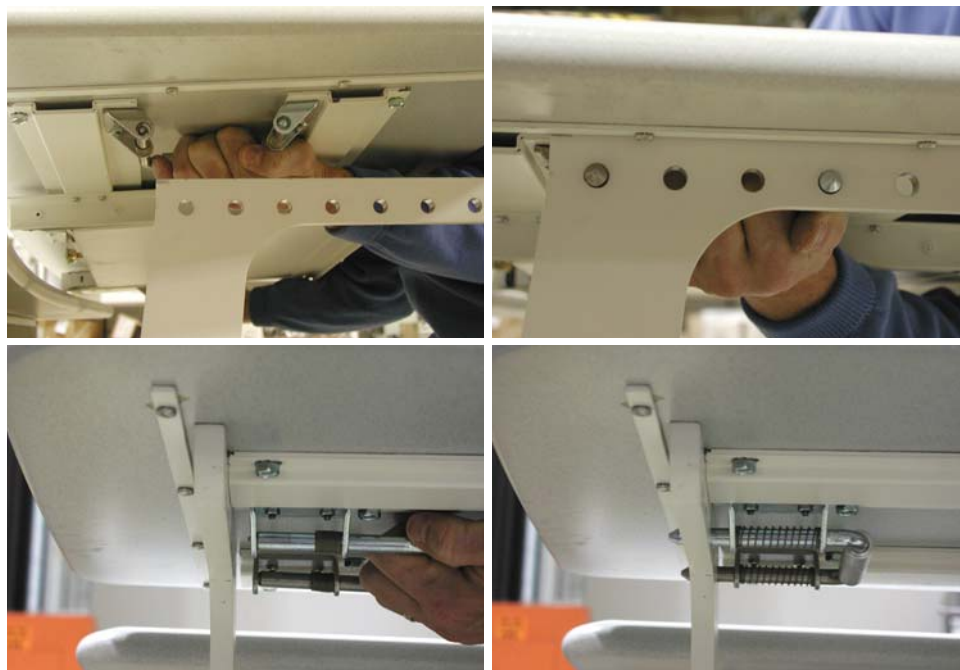


- 1.) Locate the spring loaded latches at either end of the underside of the keyboard table top.
- 2.) Pull latches inward (toward center of table top).
- 3.) Position table on brackets so that latches align with holes.
- 4.) Release latches, making certain that the pins fully engage the bracket holes (pins should protrude beyond outside edge of bracket).

9.5 Monitor and SCSI Tower Placement

- 1.) Locate and unpack the two monitors.
- 2.) Carefully place the monitors onto the monitor arms the console desktop.
- 3.) Locate and unpack the SCSI drive tower.
- 4.) Place the SCSI drive tower on the console desktop. A suggested placement is illustrated in [Figure 1-59](#).

Figure 1-59 Global Console tabletop component placement

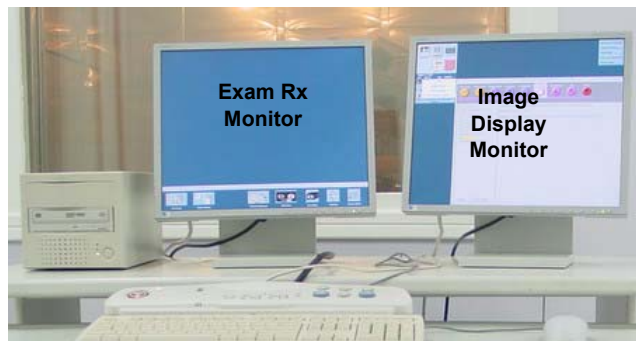

Section 10.0

Install Operator Console With FWS Table

10.1 Unpack Console

- 1.) Remove all items from the console.
- 2.) Remove all packing materials and discard.

Note: For the Operator Console shipped without skid package, go to [10.2](#).

- 3.) Place the step-board under the front edge of the skid and step on it to raise the front edge of the skid as in [Figure 1-51](#).

Figure 1-60 Step-board used to raise front edge of skid



- 4.) Remove the two front cushions from the bottom of the skid. Refer to [Figure 1-52](#)

Figure 1-61 Cushion on bottom of skid



- 5.) Seismic Brackets are shipped with system. Remove the Seismic Brackets from each side of the console. Refer to [Figure 1-53](#). Save brackets for use later if you need to mount the console

to the floor.


Figure 1-62 Removing the Seismic Brackets



- 6.) Lift up on the strap on the front of the step-board ([Figure 1-51](#)) to lower the skid. Remove the step-board.
- 7.) Ensure the console stabilizers are in line with the notched portion at the front of the skid. This will allow enough clearance to smoothly roll the console down the ramps.
- 8.) Move console to installation location.

10.2 Remove Console Covers

Figure 1-63 Front Cover Removal - Global Consoles

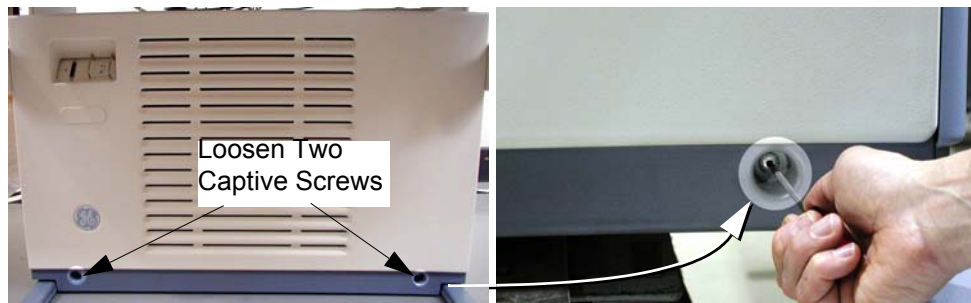

- 1.) Loosen two captive screws at bottom of console.
- 2.) Rotate bottom of cover outward and upward until it may be lifted free of the console at the top.

10.3 Replace Console Table Top With Black Table Top

- 1.) Refer to [Figure 1-64](#). The console normally arrives with the bottom hole on the console table top aligned with the fourth bolt from the bottom.

Figure 1-64 Console Table Top



- 2.) Replace the Monitor table top with the black table top:
 - a.) Remove four (4) screws with 6mm hex tool at both right and left side of the console to remove monitor table top together with the support bracket.
 - b.) Put the table top on the ground with the bracket upward. Remove six (6) screws with 6mm hex tool to remove the bracket from the table top.
 - c.) Put the black table top on the ground with table surface downward. Install the bracket on black table top with original six (6) screws.
 - d.) Install the black table top on the console with the revised order in step a).
- 3.) Install the console side covers. Side covers are left and right specific:
2315656-6: Freedom workspace left side cover for console chassis
2315656-5: Freedom workspace right side cover for console chassis
- 4.) Move the console into its final position, maintaining a 6" (15 cm) minimum distance between the console and the wall.

10.4 Install FWS table and Monitor Arms

- 1.) Assemble FWS table.
Refer to [Appendix B, on page 135](#) for details of FWS and monitor arms assembly.
- 2.) Install LCD Monitor and install monitor arms.
- 3.) Place the FWS table at one side of the console.

NOTICE FWS table should be place side by side with the console considering that the extended cable length used between FWS and Console is limited within 3 metre.

10.5 Peripherals Placement

- 1.) Place keyboard, SCIM
- 2.) Locate and unpack the SCSI drive tower.
- 3.) Place the SCSI drive tower **on the console desktop**.

Section 11.0


Seismic Mounting

Seismic kits/brackets for console and Power Distribution Unit are available shipped with the subsystem or in Shipping Collector.

11.1 Console

If site specifications require seismic mounting, use ½" bolts to mount the brackets to the floor. Refer to [Figure 1-65](#) for hole placement. The console seismic brackets are shipped with the console. They are used on the console shipping skid.


Figure 1-65 Seismic Console Mounting Hole Locations

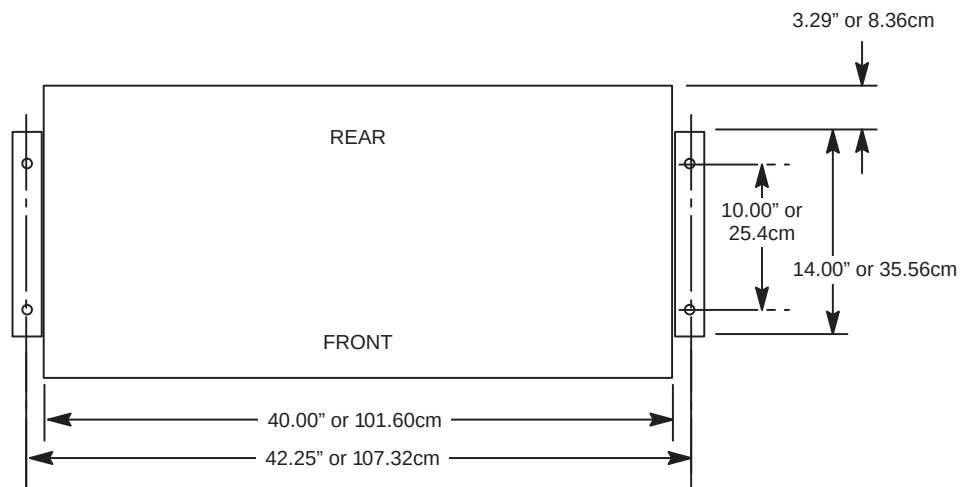

11.2 Power Distribution Unit

If site specifications require seismic mounting, use the PDU seismic brackets that were shipped with the seismic kit. Refer to [Figure 1-66](#) for hole placement.

Figure 1-66 Seismic PDU Mounting Hole Locations

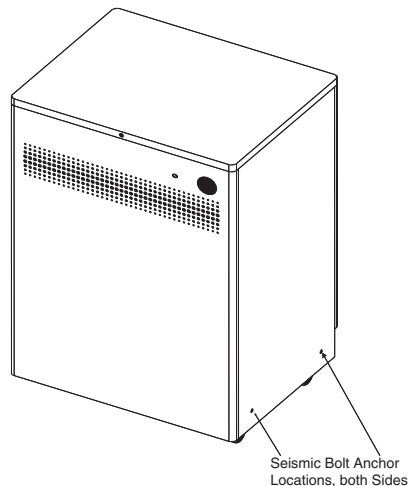

Chapter 2

Power, Ground & Interconnect Cables



NOTICE Potential for Data Loss and/or Equipment Damage

To prevent potential data loss and equipment damage, please do the following:

- When instructed, record data collected from the procedures in this chapter on GE Form e4879. For more information about this form, see [Section 8.0](#) of Chapter 4.
- Only use the Installation manual that arrives with your system for installation. Any other revisions of this manual may not exactly match your system.

Note: Use dry cleaning for electro components

Section 1.0 Introduction

Site use of conduit, floor duct, wall duct, or a raised computer floor, as well as the individual component layout determines the system cable sequence. If your site has floor or wall ducts that will interfere with placement of the table/gantry, it may be important to have the movers unload the cable boxes (8 & 9) first and run those cables while others unload the subsystems.

- Try to run the system cables after the contractor completes the contractor supplied wiring.
- All ground wires and other contractor wiring should be complete to the point of equipment placement.




NOTICE Potential for Equipment Damage

Do not store excess cable in the bottom of the PDU or Gantry.

When possible, store excess cable length in a serpentine configuration, approximately one meter long. (**Do NOT coil excess cable.**)

Figure 2-1 Excess Cable Storage Configuration

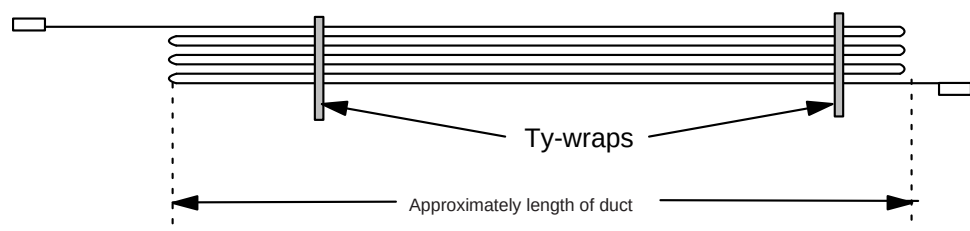

- Keep signal and control cables away from power cables and power wiring. When you lay cables in a raceway, locate the signal cables in a separate section of the raceway, or a separate conduit.
- Check all connections for tightness.
 - Use suitable tools and judgment.
 - Check all visible connections, especially ground connections.
- Check for reasonable cable routing.
 - Take into consideration necessary take-up distances for equipment maintenance, etc.
 - Try to complete as neat a job as possible.

1.1 System Component Identification

Identify all system cables by the system component designators listed in Table 2-1 on page 74. Each end of a system cable has a label, and may have a color near the connector, (refer to Table 2-2 on page 74) to indicate the component and the jack identifier of the component.

Table 2-1 System Component Identifiers

DESIGNATOR	SYSTEM COMPONENT
CT2	Gantry
CT1	Patient Table
PM	Power Distribution Unit
OC	Operator Console (Console Computer)
WL	X-Ray ON Warning Light

1.2 Cable Color Identifiers

The ends of the cables may be marked with a piece of blue, yellow, red, or orange colored tape to help with the cable installation. Table 2-2 on page 74 lists the subcomponent, and corresponding color.

Table 2-2 Cable Color Identifiers

SUBCOMPONENT	COLOR
Gantry	Blue
Table	Yellow
PDU	Red
Console Computer	Orange

Table 2-3 BrightSpeed Select System Interconnect Cables

RUN NO.	DESCRIPTION	PART NUMBER	
		LONG CABLES (KIT 2281840-4)	SHORT CABLES (KIT 2281840-5)
1	Facility MDP to Room Disconnect (A1)	cust. supplied	cust. supplied
2	Room Disconnect (A1) to PDU	cust. supplied	cust. supplied
3	Room Disconnect (A1) to System E-Off	cust. supplied	cust. supplied
4	PDU to Room Warning Light(s)	cust. supplied	cust. supplied
5	PDU to Scan Room Door Switch	cust. supplied	cust. supplied
50	HVDC Power Cable - PDU to Gantry	2343529	2343529-2
51	HVAC Power Cable - PDU to Gantry	2343530	2343530-2
52	LVAC Power Cable - PDU to Gantry	2343528	2343528-2
53	LVAC Power Cable - PDU To Operator's Console	2343531	2343531-2
54	LVAC Power Cable - Gantry to Table	n/a	n/a
55	Ground, PDU to Raceway	2371450	2371450-2
56	Ground, Raceway to Console	2371450-3	2371450-4

Table 2-3 BrightSpeed Select System Interconnect Cables (Continued)

RUN NO.	DESCRIPTION	PART NUMBER	
		LONG CABLES (KIT 2281840-4)	SHORT CABLES (KIT 2281840-5)
60	LVAC Power Cable - PDU to Optional UPS	-	-
61	LVAC Power Cable - UPS Disconnect Panel to PDU	-	-
90	LVAC Power Cable - PDU to PET	-	-
100	Signal Cable - Gantry to PDU	5120646	5120646-2
101	Signal Cable - Gantry to Console	5120645	5120645-2
102	Signal Cable (Ethernet) - Gantry to Console	2373436-2	2373436-3
103	Data Cable (Fiber Optic) - Gantry to Console	2117848-2	2117848-7
104	Signal Cable - Gantry to Table	n/a	n/a
110	Signal Cable - UPS Control to Room Disconnect (A1)	-	-
111	Signal Cable - UPS Control to UPS Disconnect Panel	-	-

Section 2.0 System Interconnect Diagram

Figure 2-2 System Interconnect Diagram


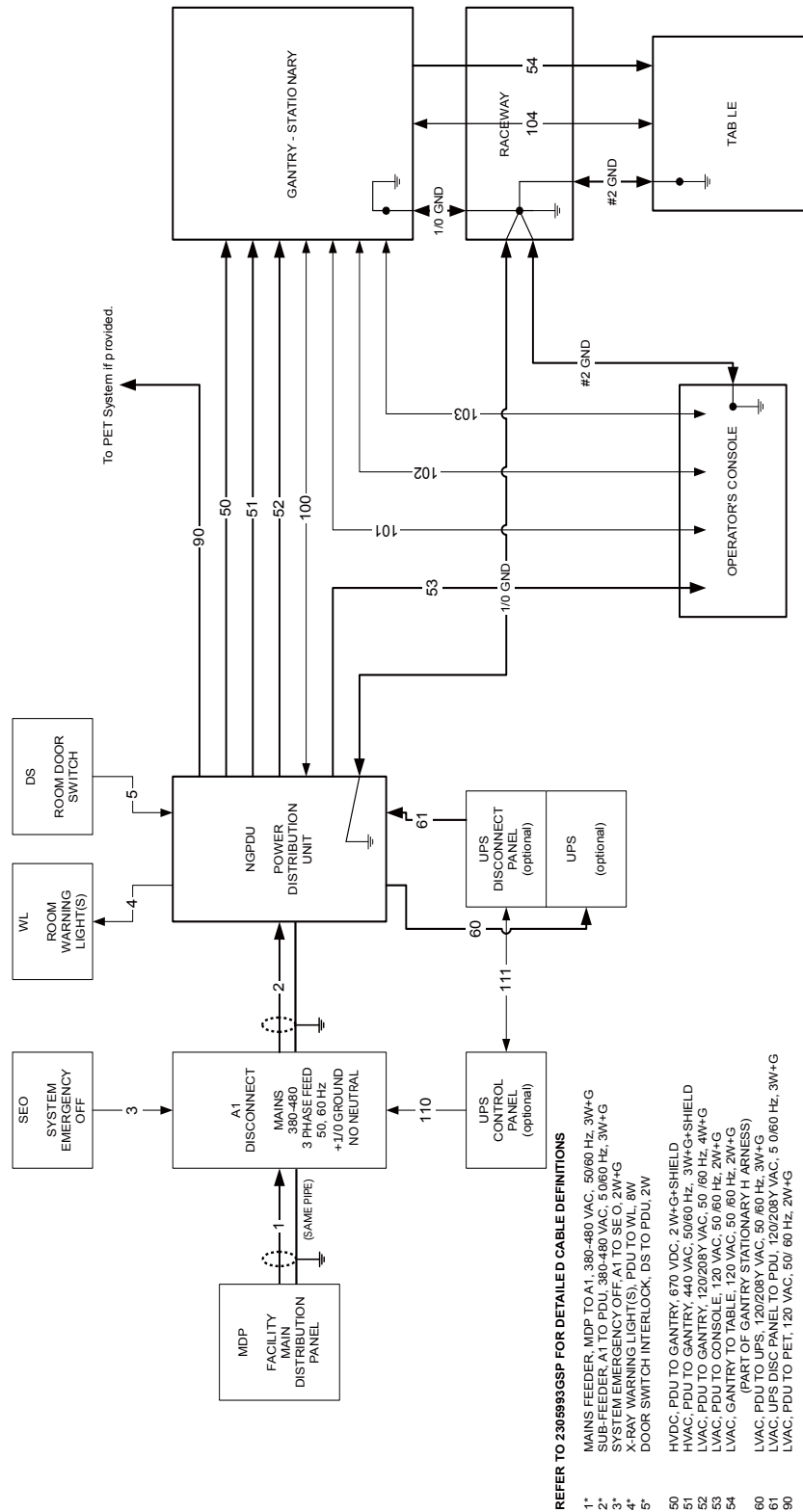


Figure 2-3 System Interconnect Diagram (with AC Outlet Box)

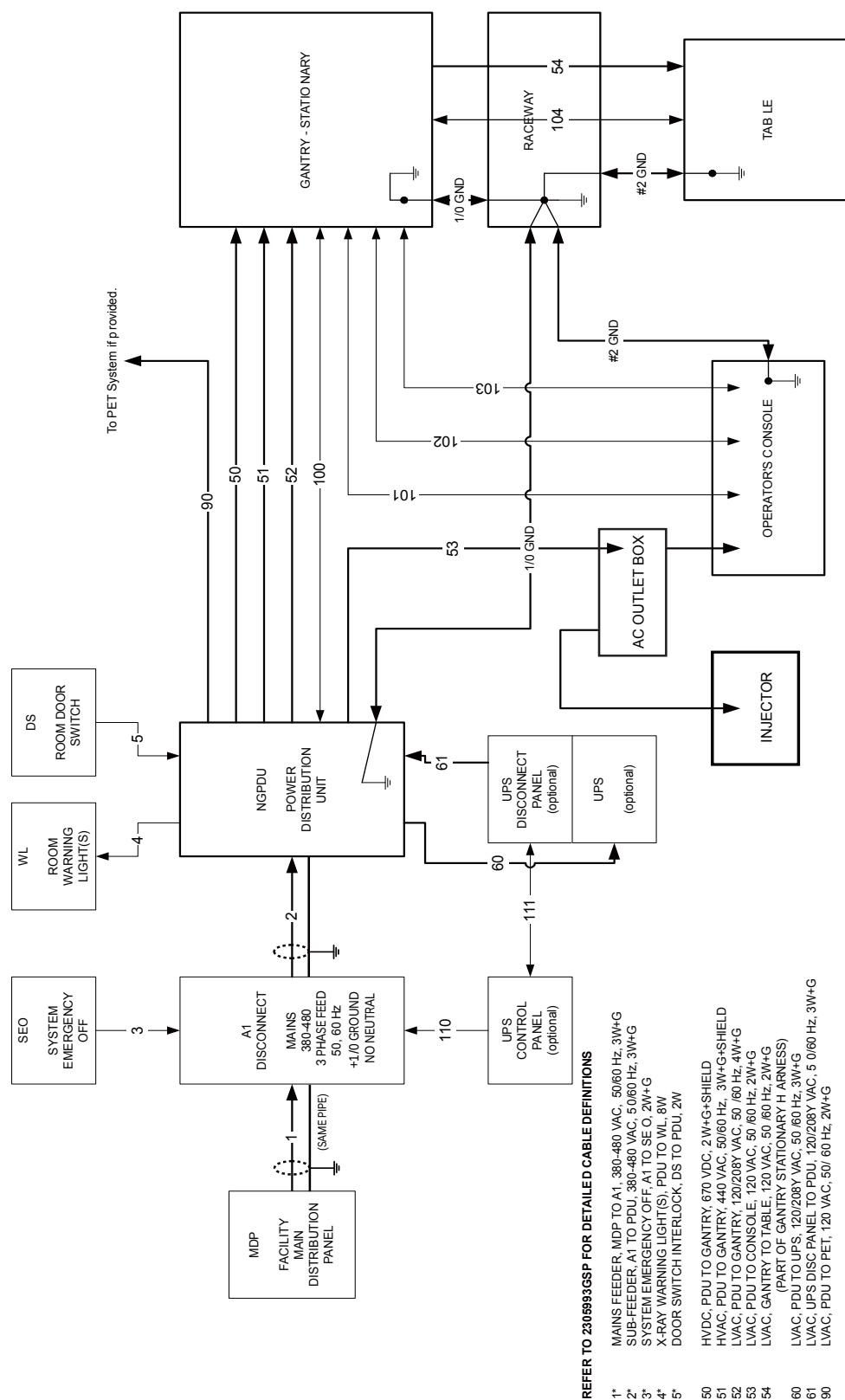

2 - Install Power

Section 3.0

Contractor Connections

Table 2-4 Contractor PDU Connections

CONNECTION OR WALL BOX	AWG #	CONNECTION FROM	CONNECTION TO PDU	INSTALLED AND CHECKED
A1	#1	Load - T1	TS-1 L1	
	#1	Load - T2	TS-1 L2	
	#1	Load - T3	TS-1 L3	
	#1/0	GND	TS-1 GND (Do NOT connect anything to neutral point.)	
WL (Warning light) See Figure 2-32 , on page 100.	#14	LV Source -1	TS6 1	
	#14	LV Source -2	TS6 2	
	#14	X-Ray ON Light -1	TS6 3	
	#14	X-Ray ON Light -2	TS6 4	
	#14	Sys-ON Light -1	TS6 5	
	#14	Sys-ON Light -2	TS6 6	
	#14	Ready Light -1	TS6 7	
	#14	Ready Light -2	TS6 8	
DS (Scan Room Door Switch)	#14	Door SW-1	TS6 9	
	#14	Door SW-2	TS6 10	

Note:
Add #2 ground
wire.

IMPORTANT: Add AWG #2 ground wire from Table frame to Table/Gantry raceway ground bar (as shown in [Figure 2-2](#)).



WARNING



WORK WITH THE ELECTRICAL CONTRACTOR TO BE SURE EXTERNAL POWER SOURCE IS TURNED OFF.

Section 4.0 Console Connections

4.1 SCIM, Keyboard, Trackball (Option) & Mouse Installation

NOTICE Use OC Long Cable Kit (P/N: 5160577) for Freedom Workspace (FWS) table.

- 1.) Route the keyboard cable under the SCIM, as shown in [Figure 2-4](#).

Figure 2-4 SCIM control with keyboard cable routed through SCIM



NOTICE



Potential for equipment damage

Never connect a mouse or keyboard with the host computer powered “ON”. Doing so can destroy components within the host computer.

- 2.) The keyboard and mouse connect to the ports on the top of the console front bulkhead. The keyboard connects to the port on the right, the mouse connects to the left port ([Figure 2-5](#)).

Figure 2-5 Global Console front bulkhead, showing keyboard and mouse connection



- 3.) Connect the SCIM cable to the SCIM as shown in [Figure 2-6](#) (note the cable routing).

Note: Make sure the SCIM connector fits snug. Some molding may need to be removed to allow the cable to fit snug.

Figure 2-6 SCIM bottom, showing cables and keyboard mounting bracket




- 4.) Select and install the proper overlay with the appropriate language for your system: (1) with Tilt or (2) without tilt.
Verify that none of the buttons get caught and stuck under the overlay. Pay close attention to the prescribed tilt button on systems with the tilt feature.
- 5.) The keyboard should attach to the SCIM using the supplied Velcro strip and fit snugly against the SCIM when finished, as shown in [Figure 2-7](#).

Figure 2-7 SCIM connected to the keyboard with the US English tilt overlay installed



4.2 Connecting the SCSI Tower

Figure 2-8 SCSI tower: Connections (left) and Optional Seismic Mounting Bracket (right)

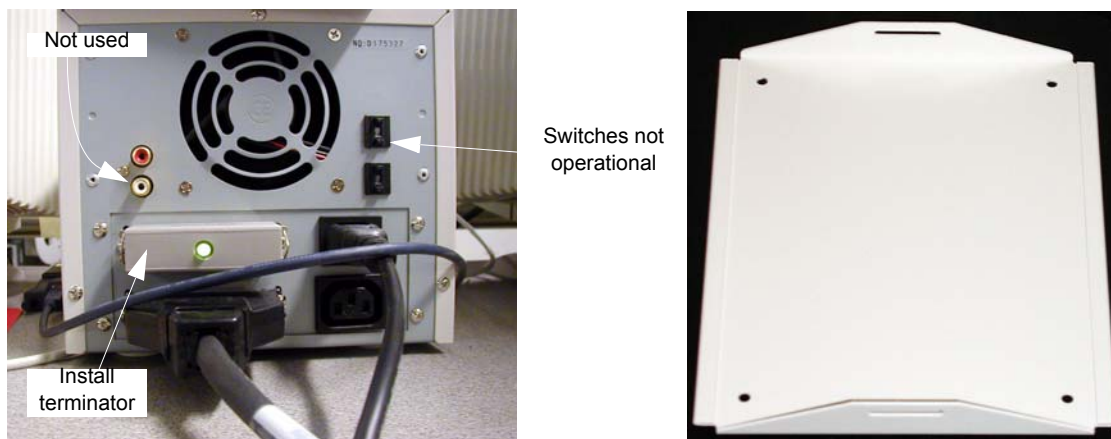

- 1.) Connect the SCSI cable to the rear of the SCSI tower.
- Note: The cable is included in the SCSI Cable Kit and is the same as the DASM cable.
- 2.) Connect the SCSI terminator to the rear of the SCSI tower.
- 3.) Connect the power cable to the rear of the SCSI tower.

4.3 Connecting the Monitor

NOTICE Use OC Long Cable Kit (P/N: 5160577) for Freedom Workspace (FWS) table.

- 1.) With the monitor and computer switched off, connect the video signal cable to the monitor's video input (see [Figure 2-9](#)).


NOTICE Equipment Damage Possible



Do not touch the video signal cable connector pins as this might bend them. When connecting the video signal cable, check the alignment of the HD15 connector. Do not force the connector in the wrong way, or the pins might bend.

- 2.) Connect the power cord to the monitor, then connect it to console power panel outlet J10 or J12.

Figure 2-9 Monitor Connections (NEC Model 1980SXi)

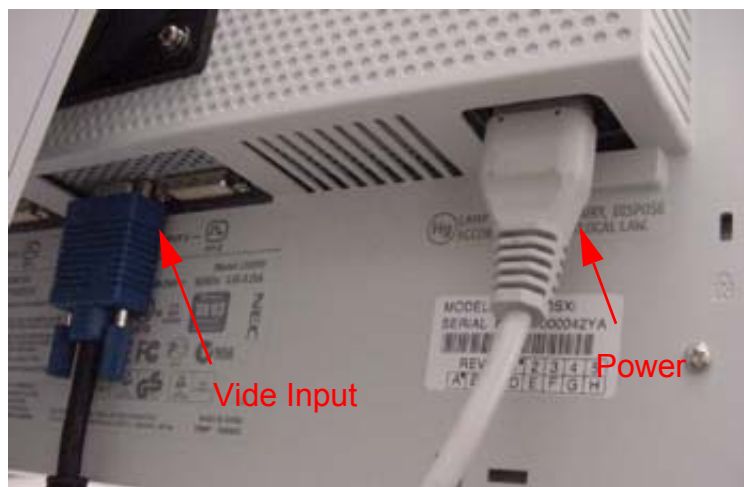

4.4 Power Panel Connections


Note: Console power is single phase power. Outlet assigned is not critical.

- 1.) Connect the console power cable to the console power panel.
- 2.) Connect console component power cords as listed in [Table 2-5](#). (“J numbers” increment from top to bottom, left to right)

Table 2-5 Power Panel Outlet Assignments - GC-GRE


J#	DEVICE	J#	DEVICE
J1		J9	Host Computer
J2		J10	Power Monitor
J3	Modem	J11	DARC
J4	IG 3	J12	Power Monitor
J5	IG 1	J13	ICOM
J6	IG 2	J14	Power Media Conv
J7		J15	Fan
J8		J16	

Figure 2-10 Power Connections

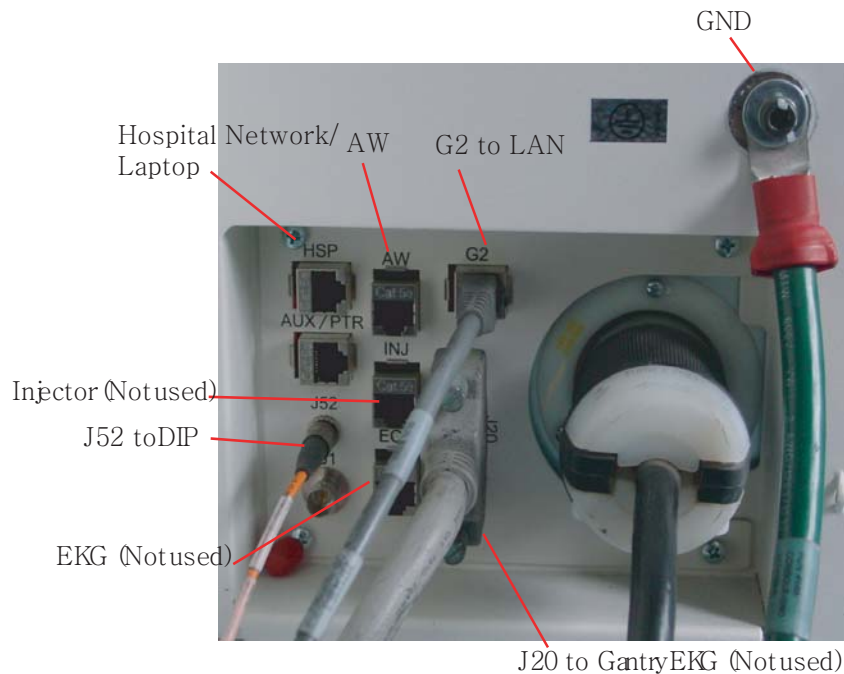

Note: If Integrated Injector is installed, AC Outlet Box must be used for power connection. Refer to [Figure 2-11](#) for cable connection.

Figure 2-11 Power Connection (with AC Outlet Box)

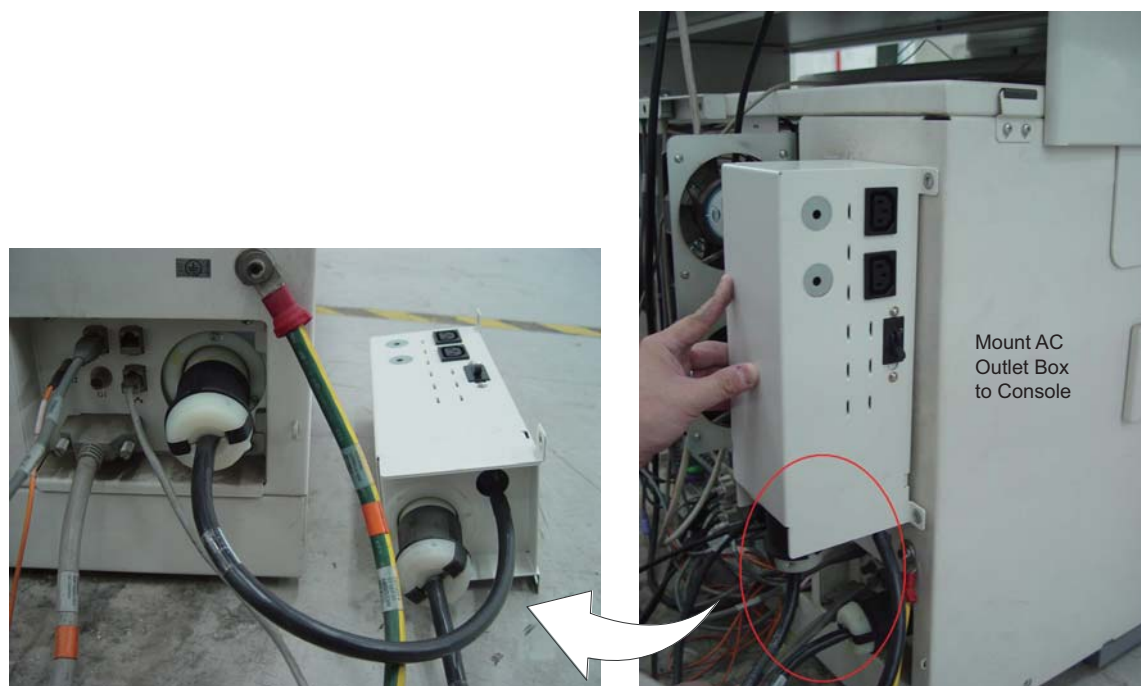


4.5 Installing a DASM (Hardware Option)

If you have a digital DASM to install, complete the following steps to connect it to the computer.


- 1.) Turn the DASM over so you are looking at the bottom.
- 2.) Set the DASM SCSI ID to 1.
- 3.) Set the DASM termination to ON.

Figure 2-12 DASM Bottom

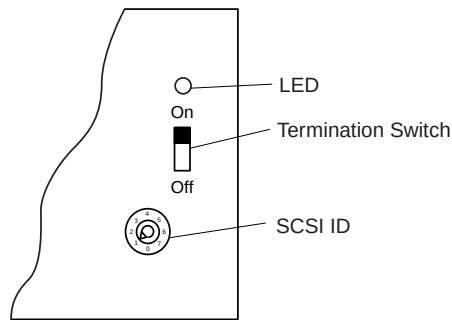

- 4.) Set the RS232/RS422 switch to the RS232 position, for a digital DASM.
- 5.) Refer to [Figure 2-13](#). Attach the SCSI cable supplied with DASM kit to the back of the DASM.

Figure 2-13 Rear view of HP xw8200, showing DASM SCSI cable connection

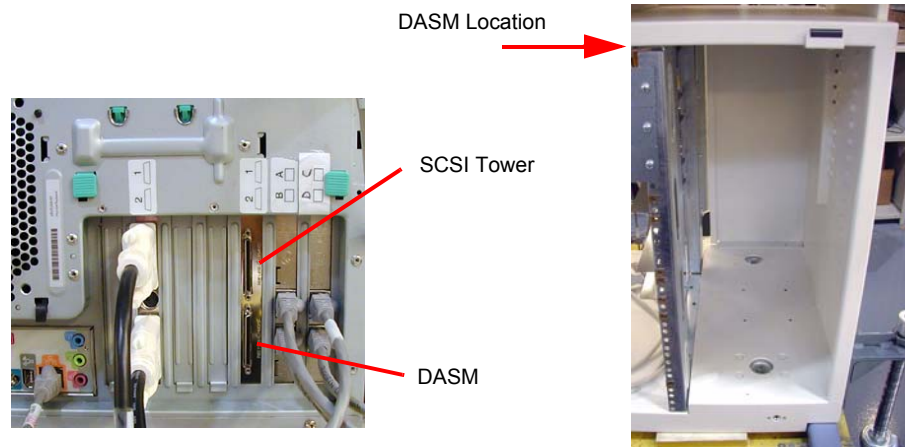

- 6.) Attach the free end of the SCSI cable to the rear of the Host Computer (see [Figure 2-13](#)).
- 7.) Plug DASM power cord into J10 (or any open slot).
- 8.) Attach remaining cable(s) from camera to DASM and install the DASM to the right of the slide-out tray.

Note: See your *Software Load Procedures* manual for instructions on configuring cameras.

4.6 Modem Option

If you have a modem to install, do it now. Place the modem on top of the console desktop. Hook up the power, phone, serial line(for GOC4/LCGOC) or USB line (for All-In-One Console) as shown in drawing which is located in the back cover of the console.

NOTICE If global Modem with Series port is used on All-In-One Console, connect the Series-USB converter cable to modem at first, then connect to USB port on the Host Computer inside of the All-In-One Console.

NOTICE If the USB Modem is used on All-In-One Console, there is no power cable for USB Modem. Only connect the USB cable to USB port on Host Computer inside of the All-In-One Console.

4.7 LAN Connections

Plug LAN cable into HSP port of rear bulkhead on console.

Figure 2-14 LAN Connection

LAN Connection



4.8 Console Interconnect

Figure 2-15 GOC4/LCGOC Interconnect


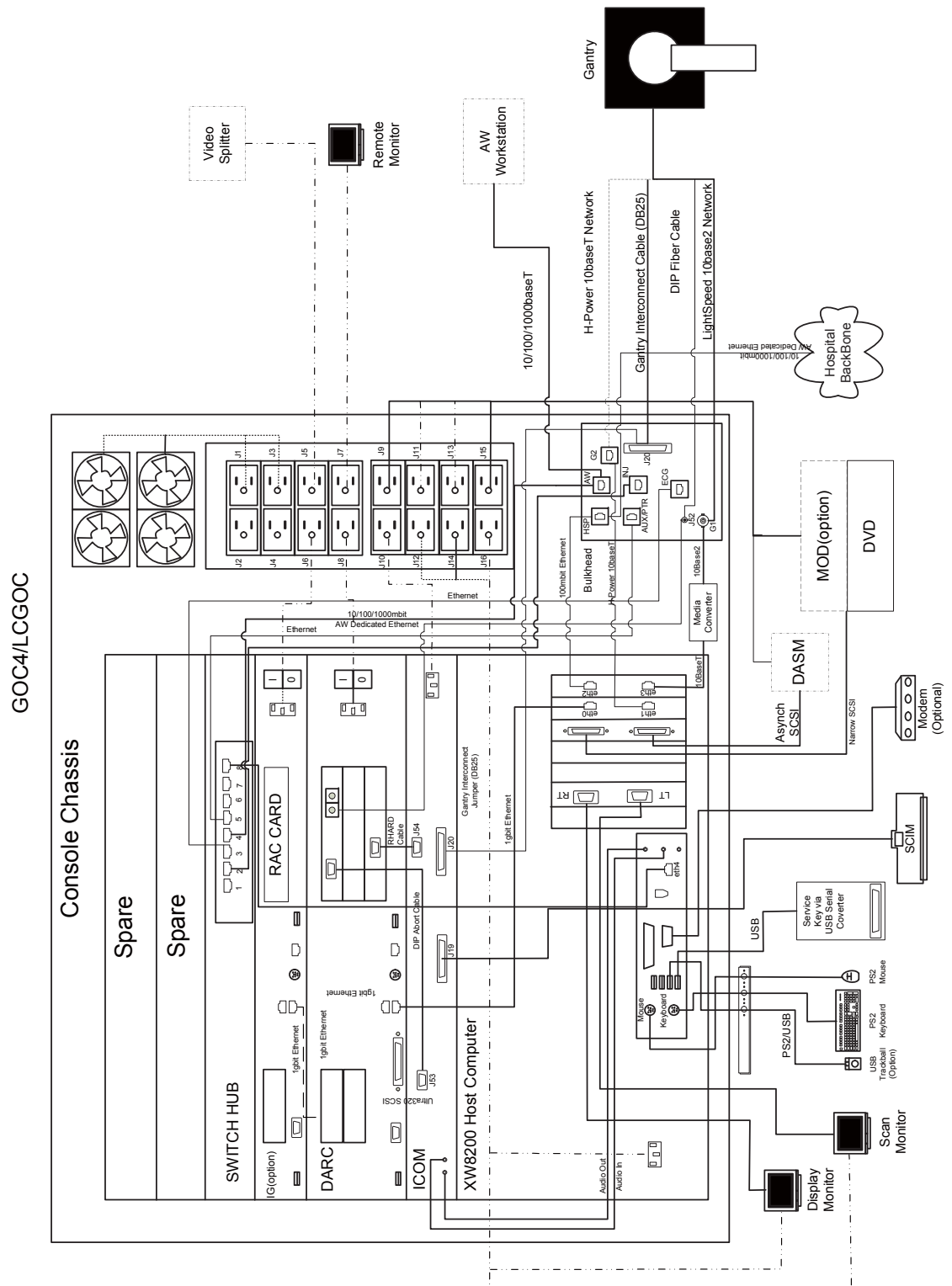


Figure 2-16 All-In-One Console Interconnect

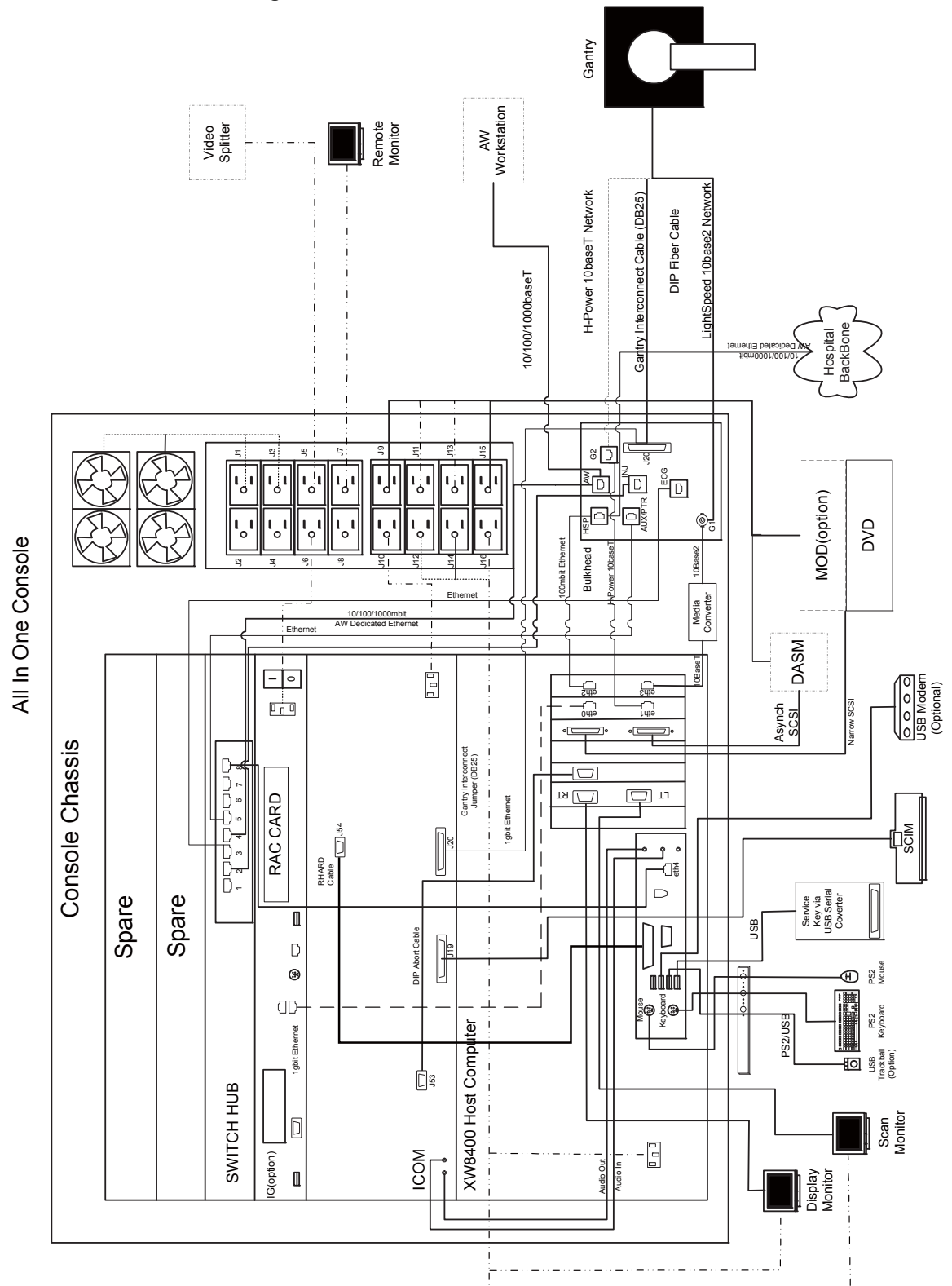

Section 5.0

Install Options

5.1 Install Optional Bar Code Reader

Follow the installation procedures in the Bar Code Reader box. When finished neatly dress all cables.

5.2 Install Optional Remote Monitor

Follow the installation procedures in the Remote Monitor box.

5.3 Install Injector Option

Follow the instructions shipped with the option. If this is a ceiling-mounted option, check that the plate is installed correctly with the holes in the correct location.

Note: For Ganty with GOB, Injector power will plug into the wall. For Gantry with IPC board, Injector power will plug into the Gantry

5.4 Install The IVY Monitor and Stand

Follow the instructions shipped with the Monitor and stand kit. Review the instructions carefully before assembling the stand and accessory basket to avoid repeated steps.

5.5 Customer Accessories (Head Holders and Extender)

Open the boxes and installed the appropriate language warning labels.

The head holders ship with shims installed to assure proper fit. Check that shims are included and a pair is installed. The holder should fit snugly. Follow the procedure in Section 4.0 of Chapter 3 to install Head Holders.

5.6 Install Service Cabinet

The service cabinet you receive may ship disassembled. Assembly takes about 1.5 hours.

- 1.) Assemble the cabinet following the instructions located in the cabinet's shipping box.
- 2.) When you complete assembly, place the cabinet in the location shown on the site print.
- 3.) Place all service materials shipped with the system in the service cabinet.

5.7 Install UPS

Follow the instructions shipped with the UPS option. The option ships with two sets of instructions, a GE set and a Powerware set.

The GE set instructs you to install the connections between the UPS and the PDU and between the UPS and the A1.

Note: A GE A1 Disconnect with UPS controls is required for this option.

The Powerware set instructs you to internally connect the batteries and do a power-up check. Refer to both manuals for additional guidance.

WARNING **LOCKOUT/TAGOUT IS REQUIRED WHEN WORKING IN THE A1 DISCONNECT.**

Section 6.0

Gantry Cable Connections

Please refer to [Figure 2-2](#) for complete system interconnect details.

Table 2-6 Gantry Cable Connections

TO	FROM	CABLE DESCRIPTION
Gantry Power Pan	PDU	HVDC
Gantry Power Pan	PDU	440VAC
Gantry Power Pan	PDU	120VAC
Gantry Power Pan	Console	Fiber - <i>Take extreme care when you install the fiber optic DAS data cable. Do not step on, kink, or sharply bend this fragile DAS cable.</i>
Gantry Power Pan	Console	LAN
TGPU (J9)	Console	Control
TGPU (J11)	PDU	Control

NOTICE



- 1.) If using a rear cable entry box (B7850RC), install it now, before routing cables to gantry.

Potential for equipment damage.

Observe correct polarity when connecting the high voltage DC power. Reversing these leads will result in serious equipment damage. The HVDC positive conductors have red insulation and are labeled “ONE.” The HVDC negative conductors have black insulation and are labeled “TWO.” Lead “ONE” must be connected to lead “ONE,” and lead “TWO” must be connected to lead “TWO.”

Observe correct phase rotation when connecting the axial motor power. Phases one, two and three should be connected top to bottom.

- 2.) Install the cables to the gantry power pan. The power pan is located on the rear of the gantry at its base. See [Figure 2-17](#), [Figure 2-18](#), and [Figure 2-19](#) for connections.

Note:

The gantry 120VAC cable may not fit under the gantry frame. Install this cable before gantry placement—or remove the power plug—to route it under the gantry.

Figure 2-17 Gantry Power Pan Connections

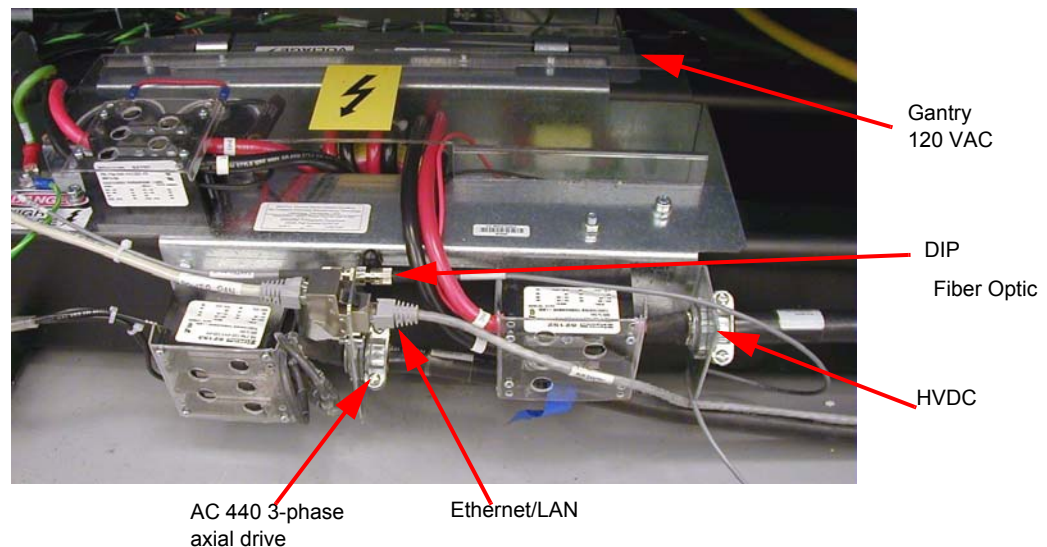


Figure 2-18 Gantry Power Pan


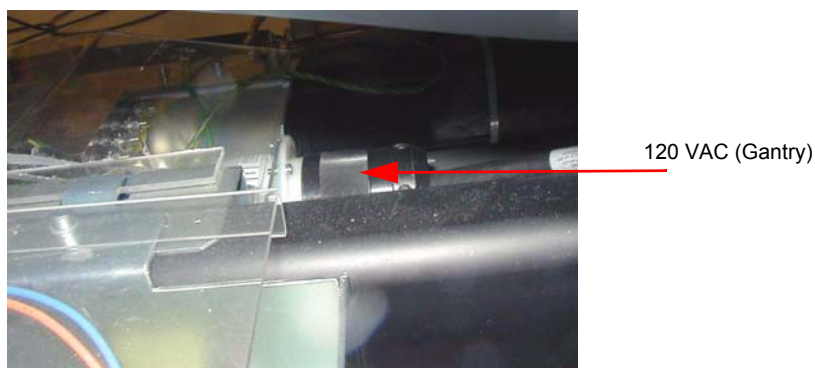


Figure 2-19 TGPU Connections

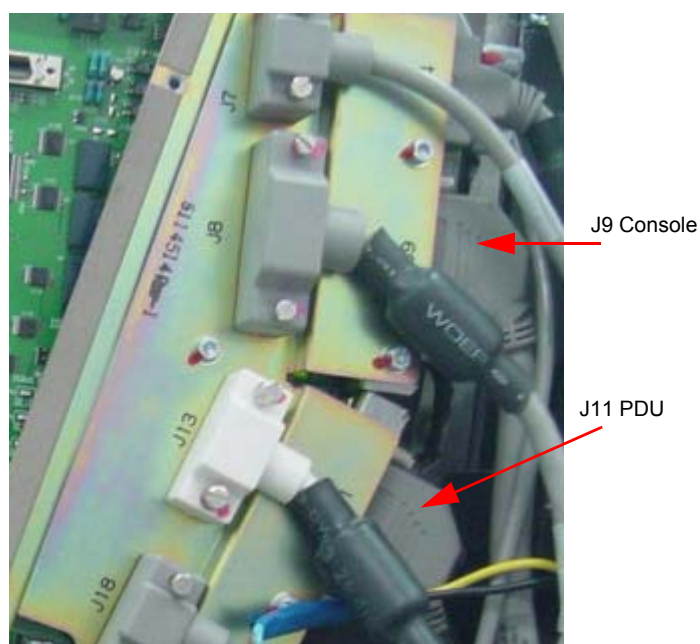


3.) Install cables to the gantry TGPU.

Note: If it's difficult to connect the cables to TGPU, follow below steps to take the TGPU outside the gantry frame (not remove the TGPU), so that it is easier to attach the cables to TGPU.

- a.) Remove the four (4) screws which secure the TGPU assembly on the gantry frame. See figure shown.

Figure 2-20 TGPU Assembly On The Gantry Frame

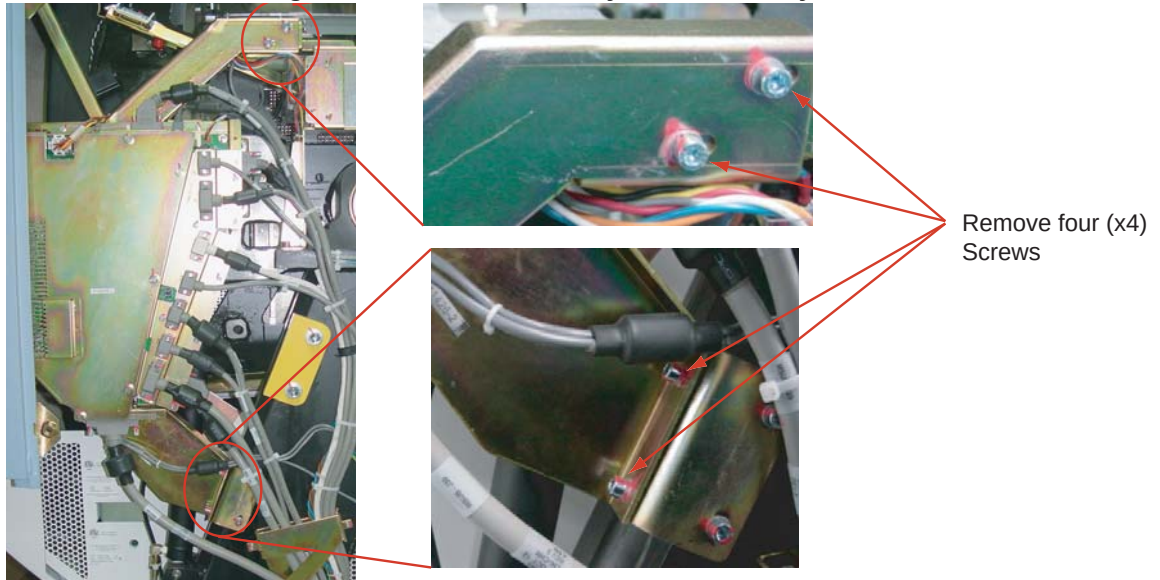

- b.) Loosen the cable clip. See Figure shown.

Figure 2-21 Cable Clip



- c.) Slightly rotate the TGPU outside, so that the cable receptacles on the TGPU is easier to access.
- d.) After connect the cables on TGPU, install the TGPU assembly with reversed order through step c) to a).

Section 7.0


Table Connections

Pull and connect the following cables:

Table 2-7 Cables Connected to Table

J#	CABLE DESCRIPTION
J1	120 VAC
J2	Signal Cable

Figure 2-22 Table Bulkhead Connections

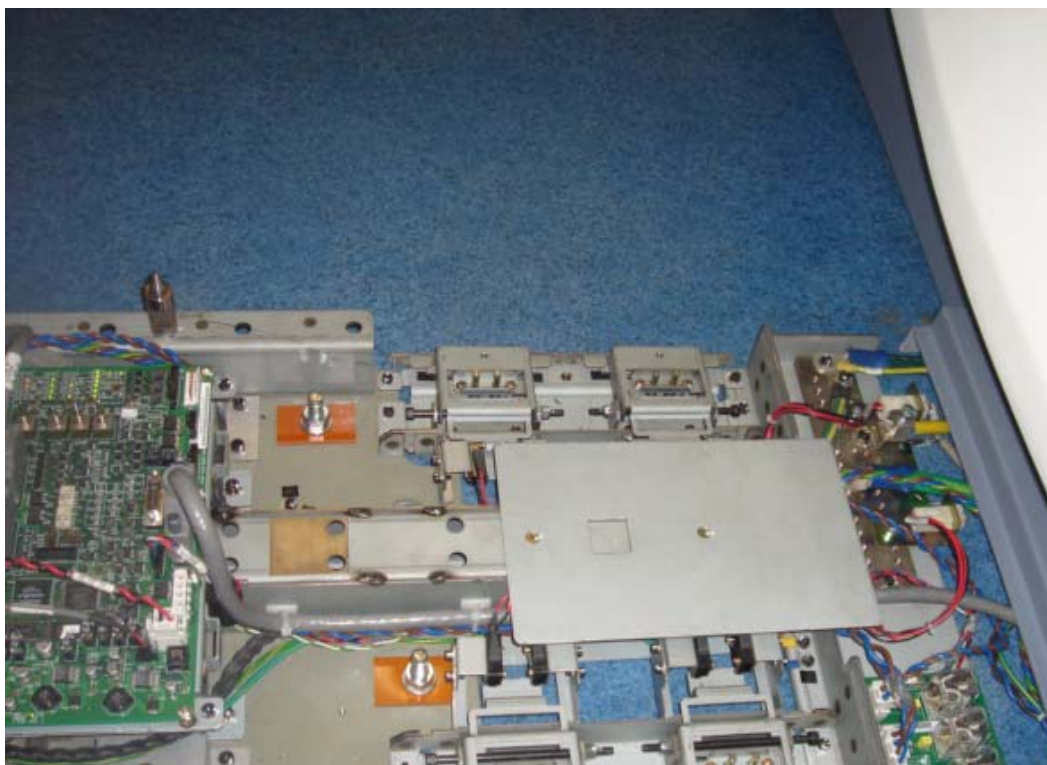

Section 8.0

PDU Cable Connections & Configuration



CAUTION

Do not work in an energized PDU. When working on the PDU, follow this simple rule: Always tag and lock out power to the PDU at the “main” disconnect. Failure to do so can result in electrocution or death.

Do apply power to the PDU until all work has been completed and all PDU covers are in their proper place.

8.1 NGPDU

As seen in [Figure 2-23](#), a number of cables must be installed throughout the PDU. Specific details on each connection can be found in the sub-sections that follow. Use [Figure 2-23](#) for reference. The PDU has been designed to have cables routed into the PDU from its behind and/or from beneath it.

Figure 2-23 PDU Cable Connections - Front


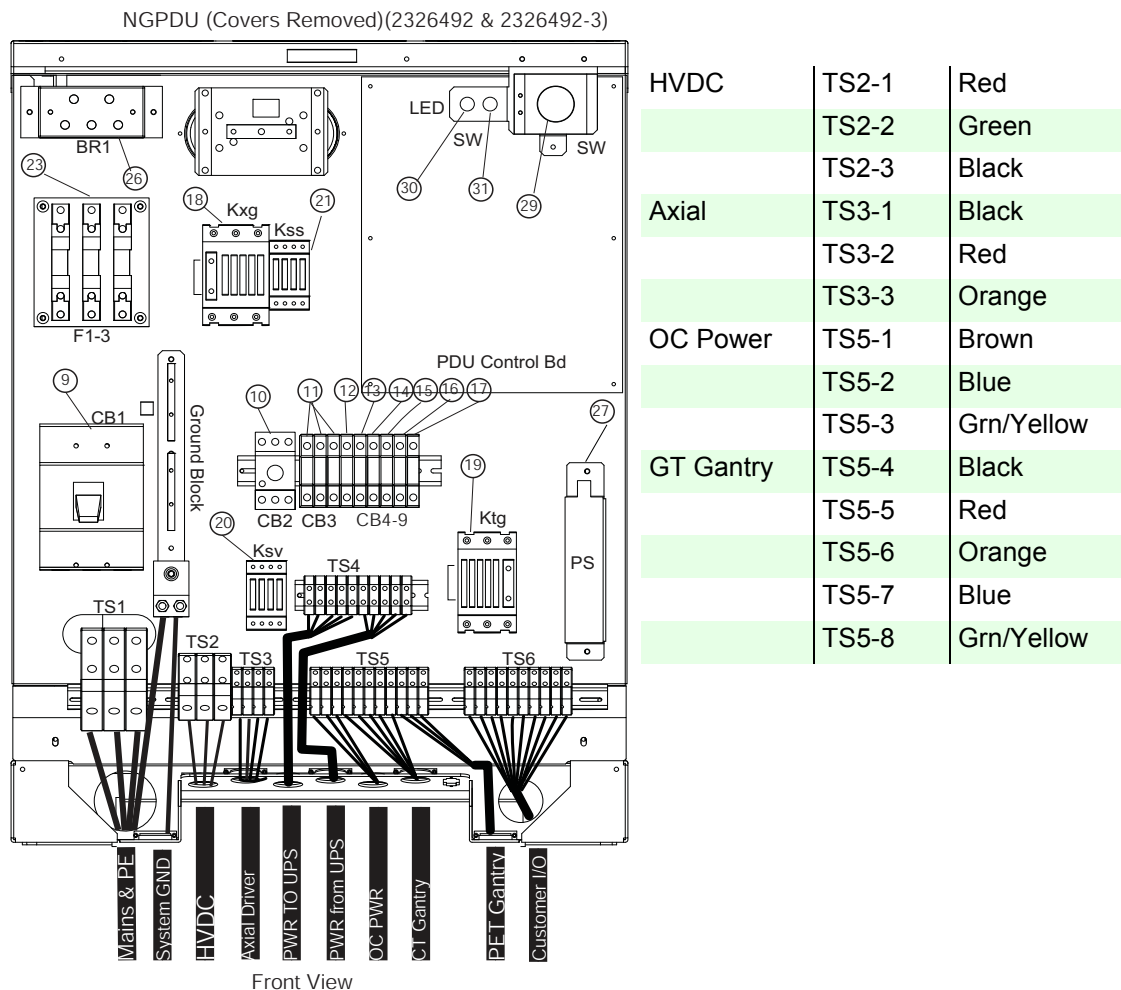



Figure 2-24 PDU Cable Connections - rear

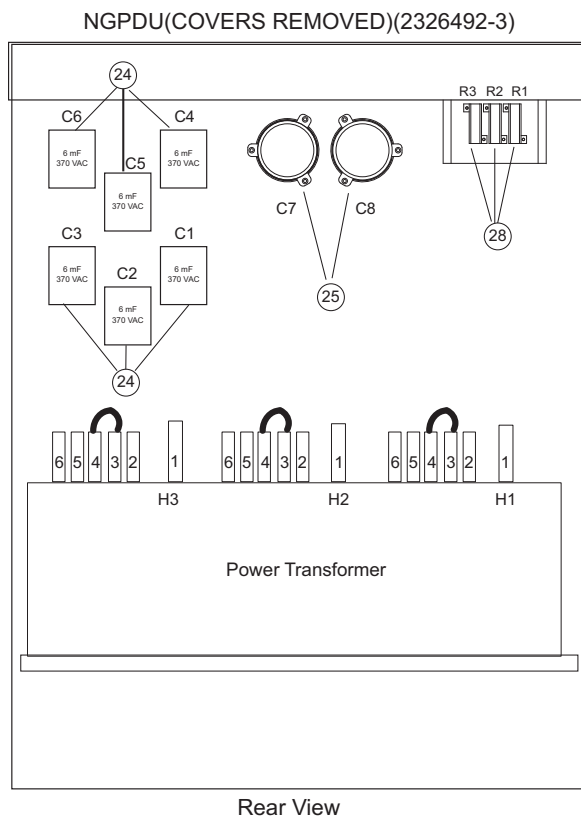

8.1.1 Panel - 380 - 480VAC Mains “TS1” Input Power Connection

- 1.) Remove the TS1 panel front cover.
- 2.) Strip the wires to fit securely on the power block.
- 3.) Observe incoming phases (L1, L2 and L3) and insert bare leads into power block.
- 4.) Insert “vault” ground into PDU “vault” ground lug.
- 5.) Tighten all fasteners securely and replace the TS1 front panel.

Figure 2-25 Input Power Panel Connections

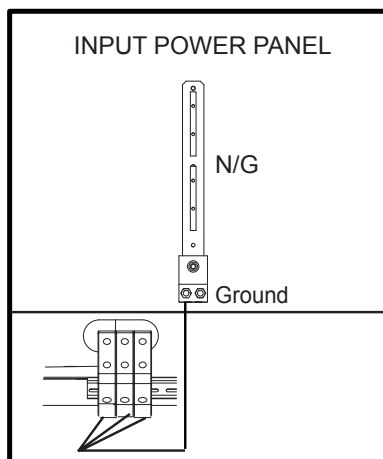


When Mains power is available to the PDU, the “TS1” power light will be illuminated.
(See [Figure 2-23.](#))

8.1.2 Panel - Circuit Breakers

Place the circuit breakers in the “off/down” position during installation, even with Mains incoming power tagged and locked out. After you have completed work on the PDU, you may return the circuit breakers to the “ON” positions.

Figure 2-26 Circuit Breaker Panel

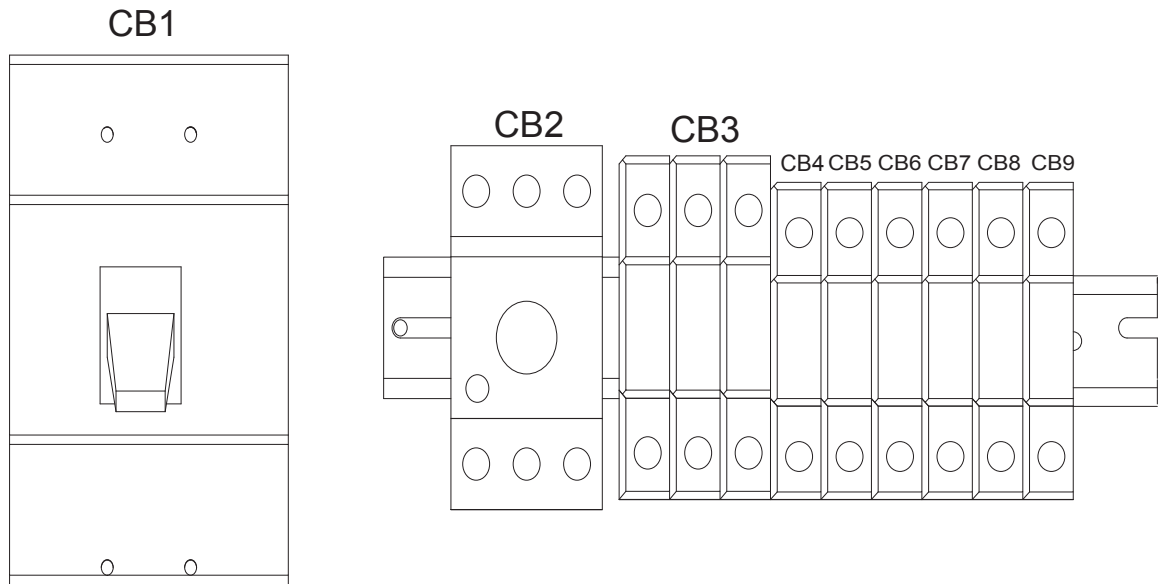

By design, when CB3 is in the “OFF” position, circuit breakers 4, 5, 6, and 7 are switched “OFF”. CB3 is essentially in series with these breakers.


Table 2-8 Panel Circuit Breaker Descriptions

CIRCUIT BREAKER	DESCRIPTION
CB3	Fully Winding Protection (Master power of CB 4, 5, 6, and 7)
CB4	CT Gantry Service Outlets
CB5	CT Gantry rotating loads
CB6	Table & CT Gantry Station Loads
CB7	Operator Console
CB8	PET Gantry
CB9	NGPDU Control Power Supply

8.1.3 Transformer (480VAC) Taps

Verify that the transformer taps are set properly. The transformer taps are set to 480 VAC operation at the factory. The taps should be set as shown in [Figure 2-27](#).

Figure 2-27 PDU Tap Positions for 480 volt operation

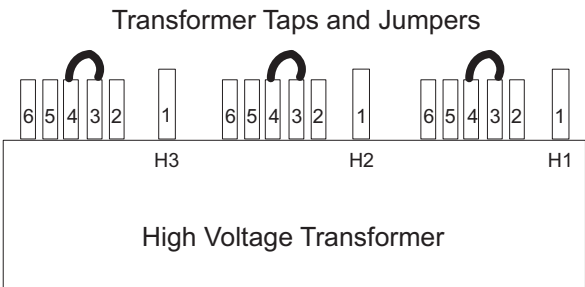


Note: Taps should be shipped as shown for 480 VAC only. For all others, you must move the taps to correct position according to [Figure 2-27](#).

8.1.4 HVDC Connection

Connect the internally shielded HVDC cable to TS2 on the standing panel. See [Figure 2-23](#) for the location of the connector and [Figure 2-28](#) for details. Observe polarities and grounds.


Figure 2-28 HVDC Connection

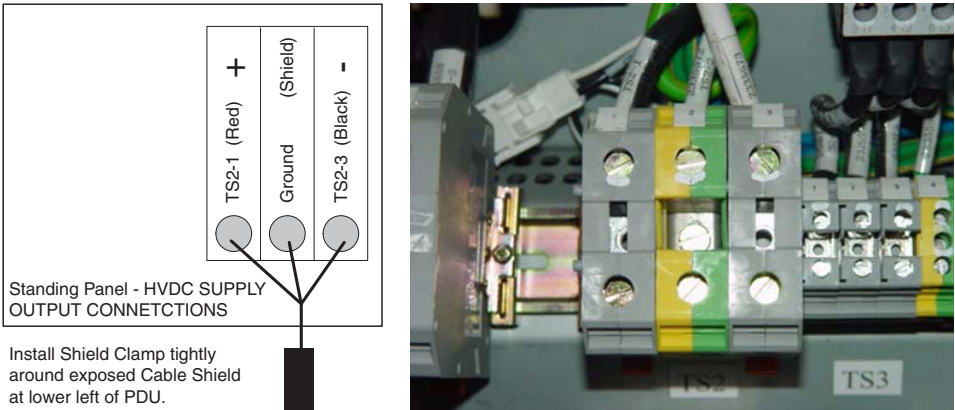

8.1.5 440V Connection

Connect the internally shielded 440V cable from the gantry to TS3 on the panel. See [Figure 2-23](#) for the location of the connector and [Figure 2-29](#) for details. Observe the labels on the cable leads for proper identification and orientation.


Figure 2-29 440VAC Connection

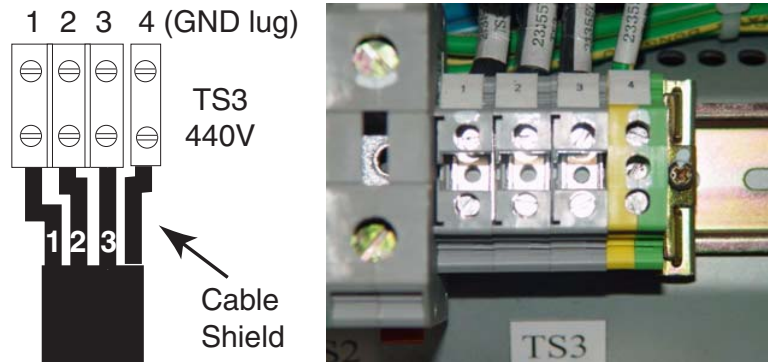

8.1.6 Gantry & Console Power Connections (120V)

Both Gantry and Console power cables come pre-terminated. Simply plug the Gantry power cable into “J4” and the Console power cable into “J5” as shown in [Figure 2-30](#).

Figure 2-30 Gantry & Console Power Connections

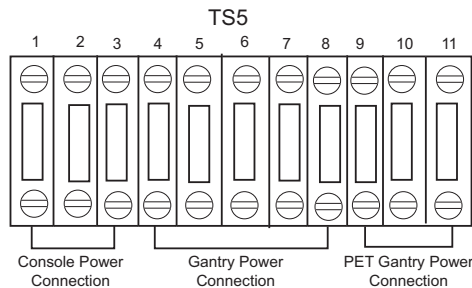

8.1.7 Console Power Cable Re-termination

Power cable re-termination should be used as a last result. Short and long cables are available.

- 1.) Carefully remove the power plug and **record the color of the wires** in [Table 2-9](#). The terminals are labelled X, W and G on the plug.

Table 2-9 Console Power Cable Termination

TERMINAL	X	W	G
Description	Hot	Neutral	Ground
Color			

- 2.) Cut the cable to desired length and dress ends.
- 3.) Re-install the power plug, according to the orientations recorded in [Table 2-10](#).
- 4.) Verify that less than 1 ohm of resistance exists between the following connections:

Table 2-10 Resistance Verification Points

FROM PDU PLUG	TO CONSOLE PLUG END	
CB1 -11 (Black) (J5 Phase X (Brown))	Phase X (Brown)	☐ Check box when complete
A3 Neutral Buss Bar (Blue) (J5 -13 Phase W)	Phase W (Blue) Neutral	☐ Check box when complete
A3 Ground Buss Bar (J5 -22 Ground Green)	Ground (Green) Ground Green screw	☐ Check box when complete

8.1.8 PDU Control Cable

The PDU control cable comes pre-terminated and should not be re-terminated in the field. Excess cable length must be stored. Simply plug the cable into “J2” on the “A4” panel. Secure it by using the fasteners integrated into cable’s connector shell.

8.1.9 System Ground Connection

Connect the ground wire (green with a yellow strip) from the Table/Gantry raceway ground bus to the system ground lug in the PDU. See [Figure 2-23, on page 94](#), and [Figure 2-31](#), below.

Figure 2-31 PDU System Ground Connection




8.1.10 Warning Light & Door Interlock Connections

8.1.10.1 Warning Light Configuration & Connection

- 1.) Warning Light is controlled by signals from the system.
- 2.) This step is site specific. The PDU by default is configured for “no” external warning light connection. If you have external warning lights, see [Figure 2-32](#) for proper connection.

Figure 2-32 Typical TS6 Warning Light & Door Interlock Connections

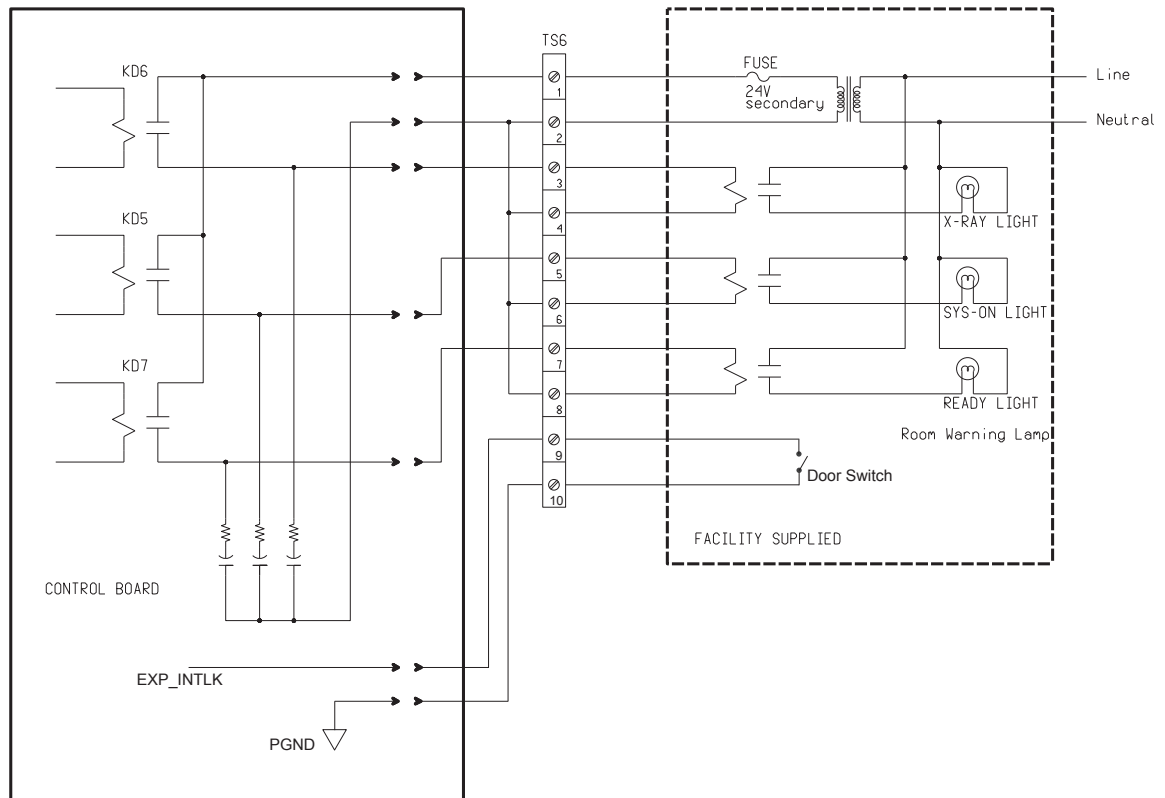


It is recommended that you use the four (4) wire method of adding a X-ray warning light to a room, as shown in [Figure 2-32](#). When using this method, you:

- Minimize EMC interference.
- Increase contact life of the relay used in the PDU.

8.1.10.2 Door Interlock Connections

Door interlocks are used to prevent X-Rays from being generated when the scan room door is open. The Door Interlock circuitry in the PDU is shipped from the factory engaged. This means the system cannot generate X-ray until disengaged. A short must exist between pins 9 & 10 for X-ray to be generated. Using a small piece of wire, short pins 1 and 2 together. See [Figure 2-33](#).


Figure 2-33 Without a Door Interlock

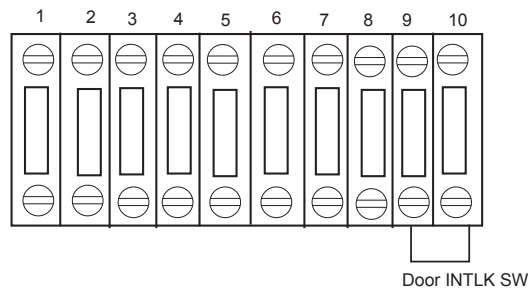

**If jumper is not in place,
exposures will not be made.
Check this jumper if you get
scan interlock errors.**

To use the system with a door interlock, wire a normally open switch between pins 1 & 2 that is attached to the interlock.

Figure 2-34 With a Door Interlock

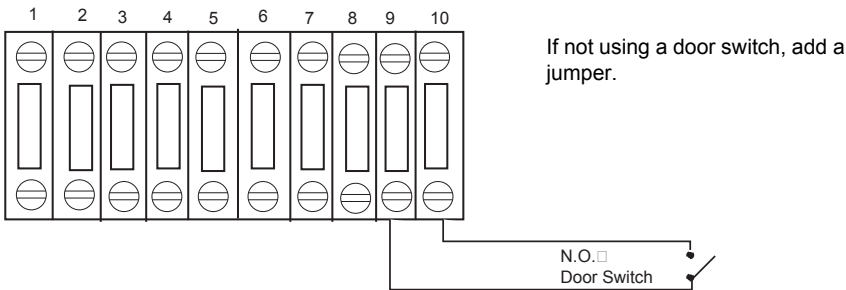


Section 9.0

System Ground Connections

As seen in [Figure 2-35](#), the Table/Gantry raceway ground bus is used to centralize all system grounding. The system ground is tied to vault ground at the PDU, through its chassis.


Figure 2-35 CT System Ground Connections

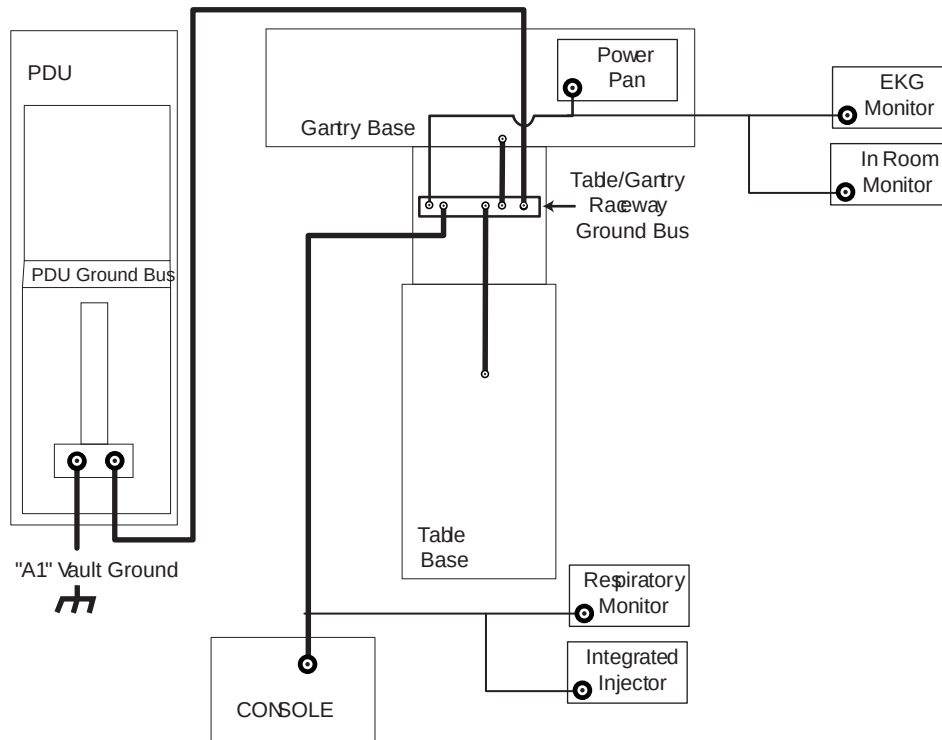

The gantry is tied to system ground at a number of points. It is important that all of these ground connections are securely made. See [Figure 2-36](#).

Check that the terminals will not be loosened, by moving or swinging the cables rather strongly by hand.

Figure 2-36 Table/Gantry Raceway Bus - Grounds

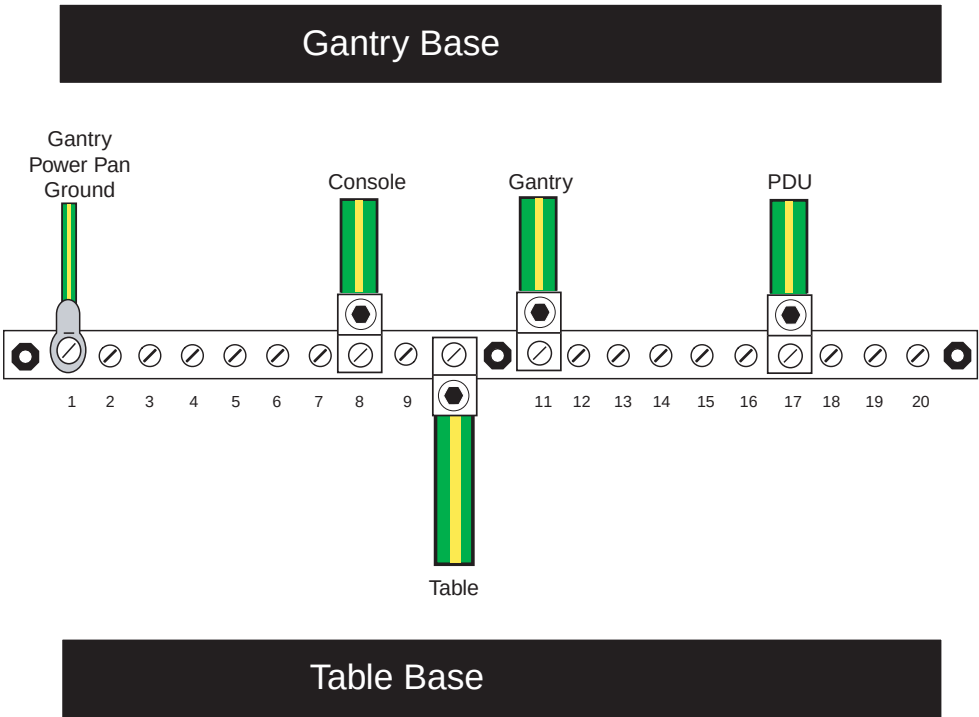

Various types and sizes of wire are used to ground the system. Please use the type and sizes specified in [Table 2-11](#), below.

Table 2-11 System Ground Connections

AWG #	CONNECTION TO	CONNECTION TO
#1/0	PDU	Power Main
#1/0	Gantry (Power Pan)	Raceway
#2	Console	Raceway
#1/0	Gantry	Raceway
#2	Table (frame)	Raceway
#1/0	PDU	Raceway

Chapter 3

System Continuity & Ground Checks



NOTICE **Potential for Data Loss and/or Equipment Damage**

To prevent potential data loss and equipment damage, please do the following:

- When instructed, record data collected from the procedures in this chapter on GE Form e4879. For more information about this form, see [Section 8.0](#) of Chapter 4.
- Only use the Installation manual that arrives with your system for installation. Any other revisions of this manual may not exactly match your system.

Section 1.0

System Continuity & Ground Checks

(Mechanical Contractor Task)

Use this section to check cable and ground connections.

1.1 Tools Required

- Digital VOM
- 10m (30 ft) of #18 wire
- 600 VAC meter leads

1.2 Procedure

Reference [Figure 3-1: Front View of NGPDU, with Covers Removed](#) on page 107 and [Figure 3-2: Gantry Power Pan](#) on page 107.



WARNING



USE AND FOLLOW LOCKOUT/TAGOUT PROCEDURES; LOCK OUT WALL POWER.

- 1.) Remove all System Power at the A1 Mains Disconnect Panel. Follow Lockout/Tagout procedures.
- 2.) Put the UPS in the Service Position.
- 3.) Remove the PDU A3 input power panel cover.
- 4.) Verify, with a voltmeter, that mains power is disconnected.
- 5.) Verify that less than 1 ohm of resistance exists between the following ground connections:

Table 3-1 Mains Connections to PDU

FROM	TO	
Wall ground connection	PDU Cabinet	☐ Check box when complete

6.) Verify that less than 1 ohm of resistance exists between the following connections:

Table 3-2 Resistance Verification Points


FROM	SIGNAL NAME (COLOR)	TO	
PDU TS2-1	+HVDC (Red)	Gantry HV Power Pan TS1-1	☐ Check box when complete
PDU TS2-2	HVDC Ground (Green/Yellow)	Gantry Power Pan Chassis	☐ Check box when complete
PDU TS2-3	-HVDC (Black)	Gantry HV Power Pan TS1-2	☐ Check box when complete
PDU Ground Bus	HVDC shield	Gantry HVDC cable shield	☐ Check box when complete
PDU TS3-1	Axial drive 440vac (Black)	Gantry HV Power Pan TS2-1	☐ Check box when complete
PDU TS3-2	Axial drive 440vac (Red)	Gantry HV Power Pan TS2-2	☐ Check box when complete
PDU TS3-3	Axial drive 440vac (Orange)	Gantry HV Power Pan TS2-3	☐ Check box when complete
PDU TS3-4	Axial drive ground (Green/Yellow)	Gantry Power Pan Chassis	☐ Check box when complete
PDU Ground Bus	Axial drive shield	Gantry 440 VAC cable shield	☐ Check box when complete
PDU TS5-1	120vac Phase A (Brown)	Console Power Plug: 	☐ Check box when complete
PDU TS5-2	120vac Neutral (Light Blue)		☐ Check box when complete
PDU TS5-3	Ground (Green/Yellow)		☐ Check box when complete
PDU TS5-4	120vac Phase A (Black)	Gantry LV Power Pan A1J1 Filter - L1	☐ Check box when complete
PDU TS5-5	120vac Phase B (Red)	Gantry LV Power Pan A1J1 Filter - L2	☐ Check box when complete
PDU TS5-6	120vac Phase C (Orange)	Gantry LV Power Pan A1J1 Filter - L3	☐ Check box when complete
PDU TS5-7	120vac Neutral (Light Blue)	Gantry LV Power Pan A1J1 Filter - N	☐ Check box when complete
PDU TS5-8	Ground (Green/Yellow)	Gantry Power Pan Chassis A1J1 Filter Ground Stud	☐ Check box when complete

Figure 3-1 Front View of NGPDU, with Covers Removed





Figure 3-2 Gantry Power Pan



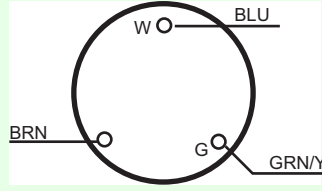

WARNING

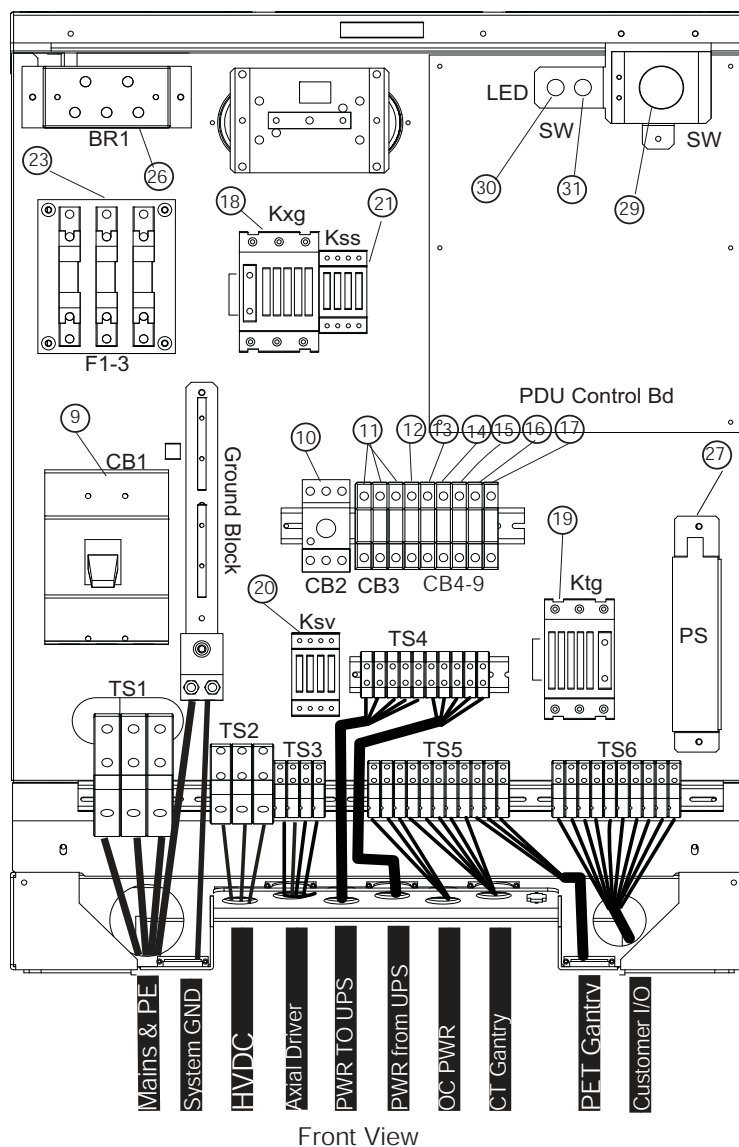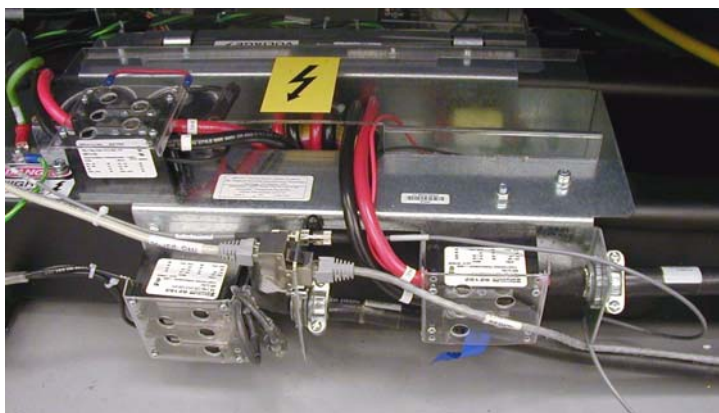

TURN OFF ALL PDU CIRCUIT BREAKERS.

- 7.) Set an ohmmeter to the lowest scale. Check between the following points for shorts to ground. Verify no continuity exists between the following points:

Table 3-3 No Continuity Verification Points

FROM PDU	TO A1 BREAKER BOX	
TS2-1 (+HVDC) (Red)	vault ground	☐ Check box when complete
TS2-3 (-HVDC) (Black)	vault ground	☐ Check box when complete
TS3-1 (440vac output) (Black)	vault ground	☐ Check box when complete
TS3-2 (440vac output) (Red)	vault ground	☐ Check box when complete
TS3-3 (440vac output) (Orange)	vault ground	☐ Check box when complete

- 8.) Leave the metal cover off the PDU A3 input power panel until you complete the checks in the next section.

Section 2.0

Site Ground Continuity Check

- 1.) Use an ohmmeter to verify the presence of **less than 1.0 ohm of resistance** between each of the following points:

Table 3-4 Resistance Verification - Site Ground

FROM	TO	
PDU Ground Bus	Vault Ground	☐ Check box when complete
PDU Ground Bus	Table/Gantry raceway ground point	☐ Check box when complete
Table/Gantry raceway ground point	Gantry Chassis	☐ Check box when complete
Table/Gantry raceway ground point	Table Chassis	☐ Check box when complete
Table/Gantry raceway ground point	Operator Console Chassis	☐ Check box when complete
All Display or Computing Options (if any)	Operator Console Chassis	☐ Check box when complete

- 2.) Install remaining covers on the:
- Gantry
 - Table
 - Raceway
 - Console
 - PDU

Section 3.0

Axial Head Holder Shim Installation

3.1 Overview

This procedure applies to following Table Types:

GT2000	All head holders
GT1700	All head holders
HP 1600	All head holders
PET Table	All head holders

3.2 Requirements

Table 3-5 Personnel Requirements

REQUIRED PERSONS	PRELIMINARY REQS	PROCEDURE	FINALIZATION
1	10 mins	15 mins	5 mins

Table 3-6 Tools and Test Equipment

ITEM	QTY	EFFECTIVITY	PART#	MANUFACTURER
Stand FE Tool Kit	1	-	-	-

Table 3-7 Replacment Parts

ITEM	QTY	EFFECTIVITY	PART#	MANUFACTURER
Shim Kit	1	-	-	-

NOTICE




Understand and Follow All General Table Safety Procedures.

3.3 Required Conditions

Check head holder for a tight fit. If the head holder fit is loose, follow this procedure.


Introduction:

- Some Axial Head Holders have a large free-play in the horizontal direction which could potentially lead to motion and therefore image artifacts.
- Installation of the 2327335 rubber shim kit can minimize this motion.

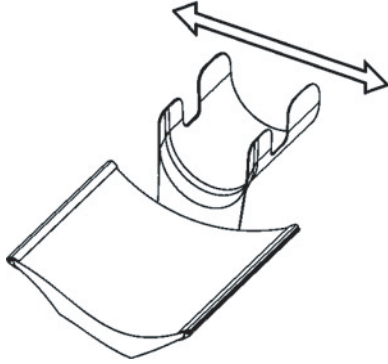


Notes before Selecting Shim Thickness:

- While selecting the best shim size, do not attach the rubber shim yet using the adhesive on the back. It is best to use a piece of tape to hold on the shim in order to see if the size is correct.
- Selecting a shim size that is too thick may result in:
 - Difficulty latching the head holder properly. The head holder must latch so that a patient is not injured.
 - Damage to the plastic latch or the plastic screws that secure it.

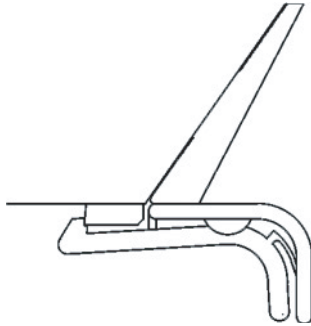


Correct - Head Holder is latched onto first step of plastic latch mechanism (The head holder does not need to be latched onto the second step)

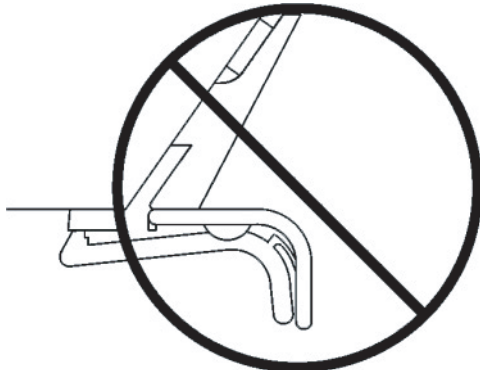


Wrong - Head Holder is NOT latched after installing shims

3.4 Procedure


- 1.) First place the two 4.0mm shims (thickest size) onto both edges of the head holder as shown (use a piece of tape to temporarily secure them)
 - The shim must be placed with the tab facing out
 - The thickness is printed on the shim

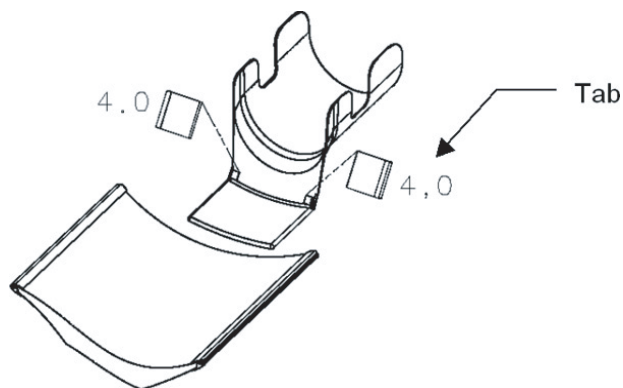

- 2.) Insert the head holder into the cradle
- 3.) Check if the head holder is latched onto the cradle at the first step of the plastic latch mechanism. (The head holder does not need to be latched onto the second step)
- 4.) Check if the head holder has a small free-play in the horizontal direction. If the rubber is too thick, repeat steps 1-4 using a thinner shim (3.5, 3.0...0.5mm) until the head holder is latched (without excessive force) and fits securely in the cradle. If the thinnest shim (0.5mm) is too tight, the tab can be cut off to reduce the thickness

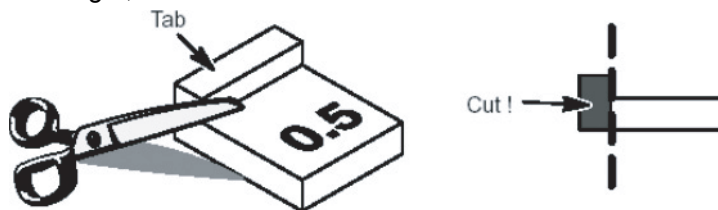

- 5.) Clean off the surfaces where the shims will mount using alcohol.
- 6.) Peel off the paper from the back of the selected shims and attach with the tabs facing out. Hold each shim with your fingers for a few seconds to attach it to the head holder.

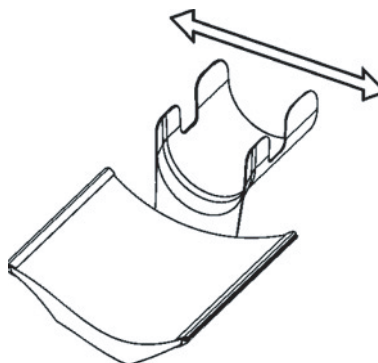

3.5 Finalization

Review latching of head holder with customer after installation.

Section 4.0

Mechanical Installation Completion Section

Complete the *Mechanical Installation* section of the GE e4879 and turn it over to your site FE. The electronic file for the checklist is found on the BrightSpeed Select Series Service Information CD-ROM 5181178-200. The electronic checklist is also located in the downloadable forms directory of the CD-ROM.

Appendix A

Removal & Installation of Covers

Section 1.0 Gantry Side Covers

1.1 Side Cover Removal

CAUTION




- 1.) Lower table to home (lowest) position.
Potential for injury if covers removed and power is left "ON".
Always remove the right side cover first, and turn OFF power at the Service Switch Panel.
- 2.) Use an 8mm Hex wrench, turn the latch $\frac{1}{4}$ turn to latch and unlatch the side cover from the front cover and the rear cover.

Figure A-1 Side Cover Latches




- 3.) Remove the right side cover by lifting it upward to release the two (2) latches, located on the top edge of the cover.

Figure A-2 Side and Top Cover Clasp

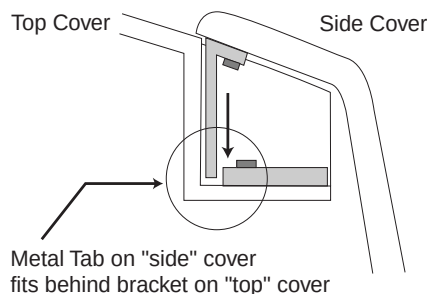

- 4.) Turn OFF the **HVDC**, **120 VAC**, and **AXIAL DRIVE** power switches on the gantry service switch panel (see [Figure A-3](#)).

Figure A-3 Gantry Service Switches

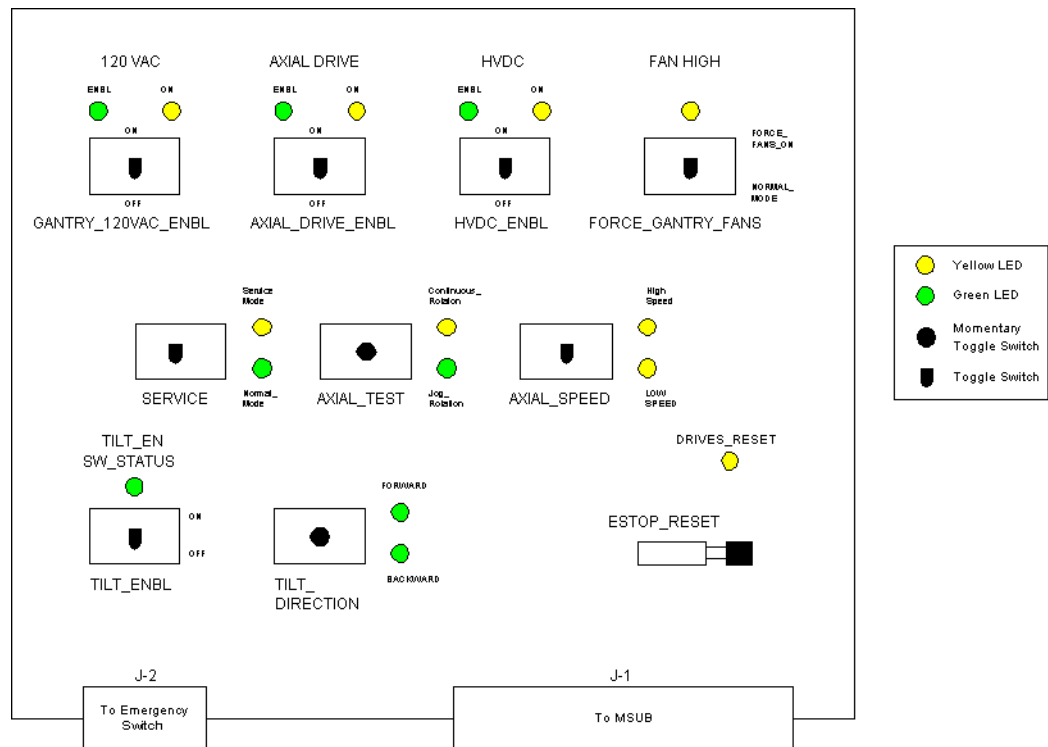

- 5.) Repeat steps 1-3 for the left side cover.

1.2 Side Cover Installation

- 1.) To install a side cover, place it over the top cover and let the two (2) side cover latches slide behind the metal tabs, located on the top cover. See [Figure A-2](#).
- 2.) Use Hex wrench to secure the side cover to front cover by turning the bolts a quarter turn. See [Figure A-1](#).

Section 2.0 Gantry Top Covers

2.1 Top Cover Removal

CAUTION



Potential for Shock.
Voltage may be present. Potential for injury if covers removed and power is left ON. Always remove the right side cover first, and turn OFF power at the service switches.

- 1.) Remove the gantry right side cover (see Gantry Side Covers Removal and Re-install).

NOTICE



Always turn OFF the HVDC before the 120 VAC. Turning OFF 120 VAC power before HVDC power can result in equipment damage.

- 2.) Turn OFF the **HVDC**, **120 VAC**, and **AXIAL DRIVE** power switches on the gantry service switch panel (see [Figure A-3](#)).
- 3.) Disconnect the power cable that connects the fan on the top cover to the gantry.
- Note: The cooling fans are now mounted on the top covers, one on each.
- 4.) Loosen the 2 Phillips screws which secure the top cover with the front and rear covers.

Figure A-4 Top Cover Mounting Screws



- 5.) Take the end of the top cover nearest to the side cover and tilt upwards.
- 6.) Slide the cover down to disengage the tab from the mounting bracket. See [Figure A-5](#).

Figure A-5 Top Cover Tabs and Bracket



- 7.) Lift the cover clear and repeat the above steps for the other cover.

2.2 Top Cover Installation

The top cover consists of two (2) pieces. Install the front and rear gantry covers, if not already installed. See [Section 3.0 on page 116](#), and [Section 4.0 on page 124](#).

- 1.) Take one of the top covers and align the tabs on the cover with its associated bracket. Lift and slide the cover into place. Position the cover to fully engage the fan interlock switch. Secure the cover using 2 Phillips screws.
- 2.) Take the other top cover and align the tabs on the cover with its associated bracket. Lift and slide the cover into place, while being sure to engage the fan interlock switch. Secure the cover using 2 Phillips screws.
- 3.) Connect the cable from the fans to the gantry.

NOTICE



Always turn OFF the HVDC before the 120 VAC. Turning OFF 120 VAC power before HVDC power can result in equipment damage.

- 4.) Turn on the three (3) power switches.
- 5.) Ensure fans work properly.
- 6.) Re-install the gantry side covers.

Section 3.0

Gantry Front Cover

NOTICE Potential for cover damage.



Front and rear cover removal and installation can be safely accomplished by one (1) person using the dollies provided with the system. Failure to use these dollies will significantly increase the likelihood of damage to the covers. Do not lean covers against walls.

3.1 Front Cover Dolly Setup

DANGER



DO NOT USE DOLLIES ON UNEVEN SURFACES SUCH AS STEPS OR ELEVATOR THRESHOLDS. THE DOLLIES ARE DESIGNED TO BE USED ON FLAT LEVEL FLOORS WITHIN THE SCANNING SUITE ONLY. MISUSE CAN RESULT IN PERSONAL INJURY OR DAMAGE TO COVERS OR OTHER FACILITY ITEMS.

CAUTION



Rotating arms on the stand are supposed to be stiff. If they fall freely, tighten the tensioning nuts. Loose rotating arms will reduce the stability of the dollies when supporting the front cover. Do not lubricate.

- 1.) Arrange Dolly sections for assembly. The base and post can be assembled only one way.
Refer to [Figure A-6](#) and [Figure A-7](#).
 - The base uses two (2) palm screws to clamp the four (4) legs in the open or usage mode.
 - The base also uses the same palm screws to prevent the legs from falling in storage mode.
 - The top post can be inserted in either base and is keyed for proper engagement.
 - The top post locking pin prevents the sections from separating during usage.

Figure A-6 Front Cover Dolly in Storage Mode


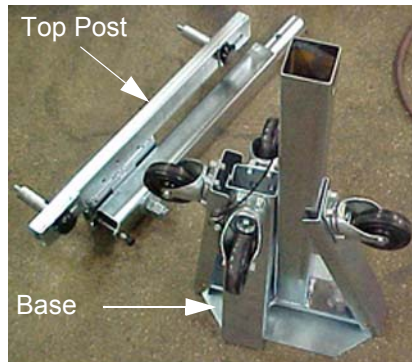


Figure A-7 Front Cover Dolly Base Assembly

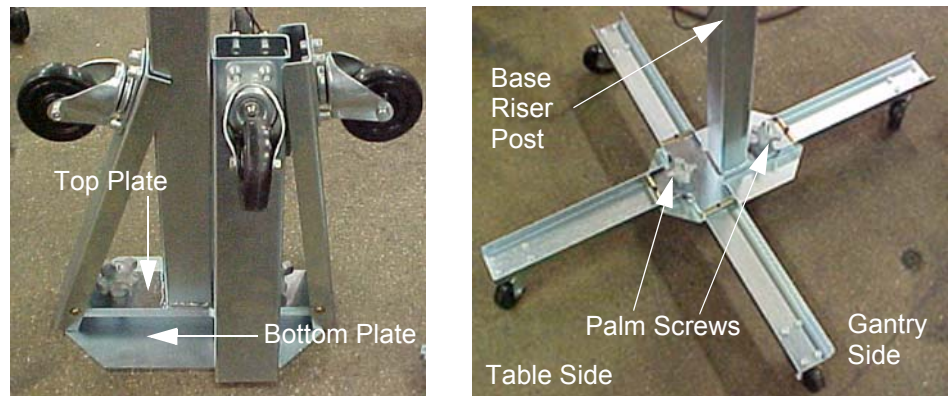

- 2.) Unfold the base legs by loosening both palm screws to the top of their travel.
- 3.) Carefully unfold the legs so that the castors touch the floor.
- 4.) Tighten the palm screws to clamp the legs between the base top and bottom plates.

Note: Lifting the base by the riser post while leaving the castors on the floor will ease palm screw tightening. Reference [Figure A-7](#).

WARNING



ENSURE BOTH PALM SCREWS ARE TIGHTENED SECURELY AND THE LEGS ARE CLAMPED TIGHTLY BETWEEN THE BASE TOP AND BOTTOM PLATES. FAILURE TO DO SO WILL RESULT IN INSTABILITY DURING FRONT COVER HANDLING.

- 5.) Insert top post into the base riser post. Align the key for complete engagement.
- 6.) Insert top post locking pin to secure both top and bottom sections.
- 7.) Reverse above steps to disassemble.


Note: For base storage only one (1) palm screw needs to be tightened. This will engage the bottom base plate and the leg ends preventing the legs from unfolding during transport and storage.

3.2 Removal

- 1.) Position the table at its lowest position.
- 2.) Remove gantry side and top covers, if you have not already done so. See [Section 1.0 on page 113](#). Make sure that the three (3) power switches have been turned off. See [Figure A-3](#).


- 3.) Assemble the front cover dolly.
 - a.) Tighten the two (2) shoulder bolts to the gantry securely. This will make cover installation easier. See [Figure A-8](#).

Figure A-8 Front Side Dolly

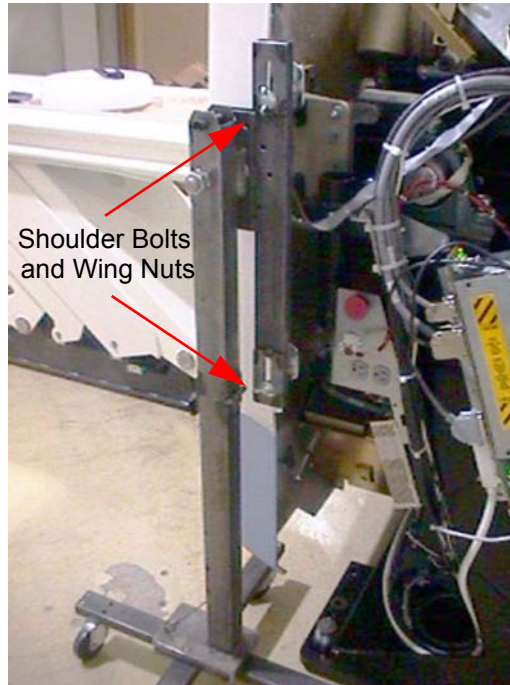

- b.) Attach side dolly to the shoulder bolts and secure assembly with two (2) wing nuts.
 - c.) Repeat steps a and b to assemble the other side dolly.
- 4.) Detach front cover J3 and J2 and front cover BKHD J1 cables.
- 5.) Remove front cover
 - a.) Disengage upper and lower cantrell brackets on both sides of the cover.
 - 1.) Using steady but firm pressure, lift each of the lower cantrell brackets from their associated retainers. See [Figure A-9](#).


Figure A-9 Releasing cover brackets

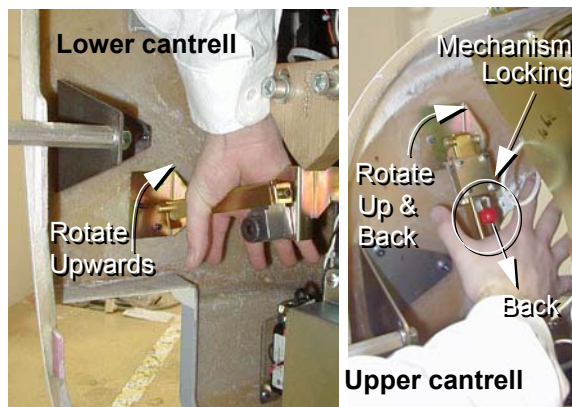

- 2.) Disengage the locking mechanism on the upper cantrell brackets by using your thumb to slide the trigger (red lever) back. This will release the locking mechanism and allow the cantrell to be rotated upwards with steady and firm pressure.

- b.) Disengage the rubber retaining straps on both sides. You may find it helpful to lift “up” on the cover to align the stud while attaching the rubber retaining straps.
- c.) Also lift and rotate cover locking arm to unlocked position.

Figure A-10 Rubber retaining straps and Cover Locking mechanism

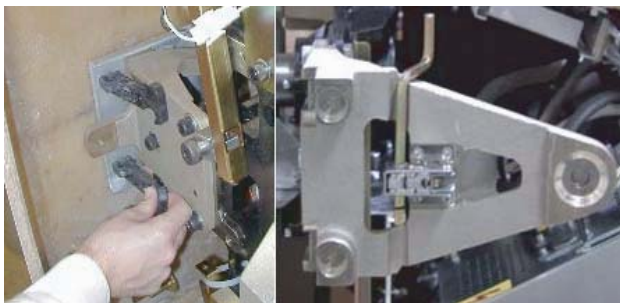

NOTE: The rubber retaining straps at the right side of the gantry are different from the current HP60 system.




- 6.) Rotate front cover away from gantry.
 - a.) Move front cover away from gantry giving enough space (about 5 feet) between front cover and gantry.
 - b.) Pull the locking pin and rotate front cover away from gantry. Place locking pin in one of the side dolly perforations. See [Figure A-11](#).

Figure A-11 releasing Front Cover Dolly Hinge



Figure A-12 Front Cover Removal Sequence

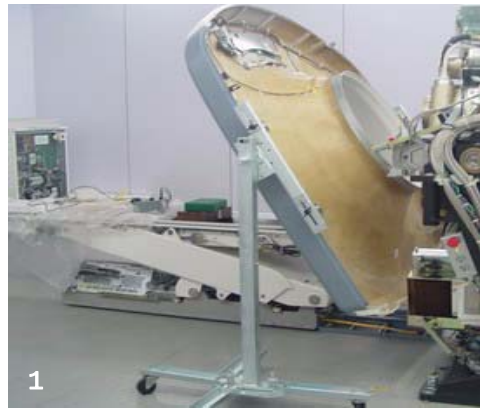

1. Upper Left>Rotate cover to make room for step 2.
2. Upper Right>Rotate cover to clear the table. Roll the cover to foot end of table
3. Lower Left> Rotate the cover upside down to provide clear work area.
Note: Remove Gantry Display and Control Panel in position 2.

- 7.) Rotate the cover horizontally and move it back and over the table to a safe location. Once in a safe location, you may over-rotate the cover full vertically but upside down.
- 8.) Remove the gantry display and one (1) of the cover's control assemblies, and place them into the service positions.
 - a.) Remove the gantry display and place it into its service position.
 - * The gantry display is held in place with (6) thumb screws. Use a flat-blade screwdriver to remove the Display. Reference [Figure A-13](#).

Figure A-13 Gantry Display Removal




- * Place the Display in the bracket on the right side of the gantry. See [Figure A-14](#).

Figure A-14 Gantry Display Service Mounting Location



- * Disconnect the cabling at the right rear gantry cover. Only (1) cable will connect to the Gantry Display. Connect the cable taken from the rear cover to the display.
- b.) Remove a gantry control and place it into its service position.
- 1.) Loosen the screws that mounts the gantry control. Keep one hand on the control panel at all times to prevent it from dropping to the floor.

Figure A-15 Gantry control panel removal

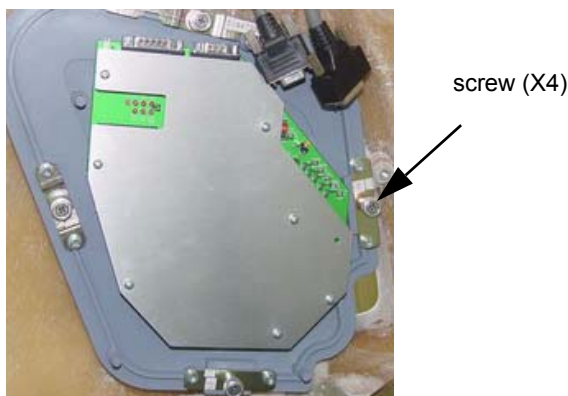

- 2.) Align the ball studs with their associated receivers and snap into place.

Figure A-16 Control panel service position



- Note:
- 3.) Connect cable to terminator located on the cantrell arm. Reference [Figure A-17](#).
There are 3 cables, each of which is unique. The ribbon cable is not used in the Service configuration. The other 2 cables will only fit in the terminator or the control panel, not both.

Figure A-17 Gantry Service Mode Cable Terminator



3.3 Installation

- 1.) Remove the gantry display and control assembly from their service positions and re-attach them to the gantry cover.
 - a.) Disconnect cables from Display and Gantry Control Panels.
 - b.) Install Gantry Display in front cover. Secure the 6 thumb screws. With a flat-blade screw driver gently tighten past finger tight.
 - c.) Install the gantry control panel making sure the screws are secure within the receivers.
 - d.) Re-attach cables.
- 2.) Rotate gantry back to its vertical position.
Potential for front cover damage.

NOTICE



When you put rotate the gantry back to its vertical position, make sure not to scratch the front cover with the edge of the table cradle.

- 3.) Attach the front cover.
 - a.) Align the studs on both sides of the front cover with each associated receiver. Receiver is located on the gantry frame.

Figure A-18 Cover stud and Mounting bracket receiver



- b.) Insert the stud on one side into its associated receiver and attach the rubber retaining straps. Then insert the stud on the other side into its associated receiver and attach its rubber retaining straps.
You may find it helpful to lift “up” on the cover to align the stud while attaching the rubber retaining straps.

- 4.) Re-attach upper and lower cantrell brackets on both sides.
 - a.) Remove upper Cantrell brackets from service position and rotate them into position over their associated retaining pins. See [Figure A-19](#).

Figure A-19 Service position of upper and lower cantrell brackets


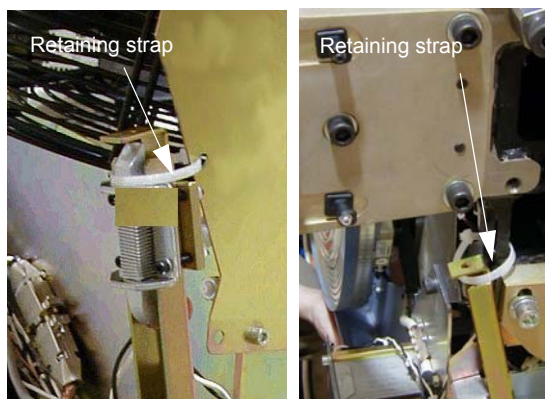


Figure A-20 Cover retaining pins (top and bottom)



Press down firmly on the bracket and snap it into place. The locking mechanism on each upper bracket should lock the bracket securely into place. Do this on both sides. See [Figure A-21](#).

Figure A-21 Locking the cover brackets into place

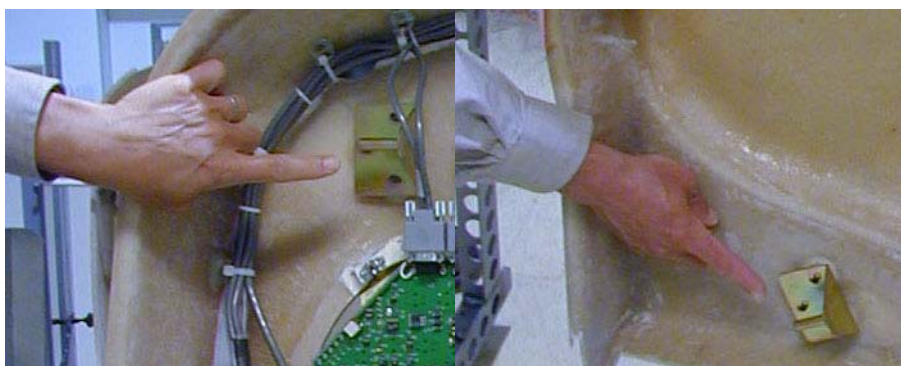

- b.) Remove lower cantrell brackets from service position (see [Figure A-19](#)), and rotate them into position over their associated retaining pins. Press down firmly on the bracket and

snap it into place. See [Figure A-21](#).

Note: Mis-adjustment of the cantrell brackets can cause misalignment of the top and side covers. The upper and lower cantrell brackets do not require adjustment during normal use.

- 5.) Remove dolly, disassemble and store safely away for later use.
- 6.) Re-attach cables to cover.


Section 4.0

Gantry Rear Cover

4.1 Removal

- 1.) Assemble the rear cover dolly.
 - a.) Tighten the two (2) shoulder bolts to the rear cover.

Figure A-22 One side of the Rear cover dolly

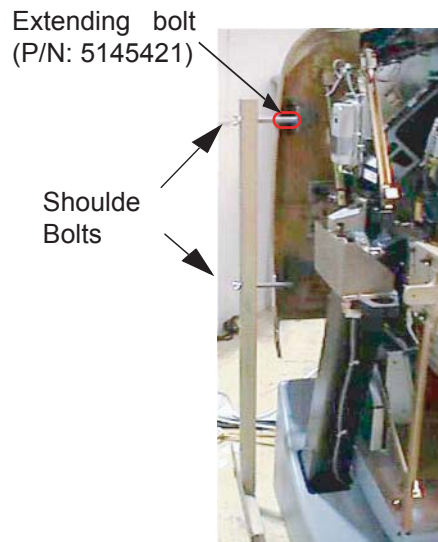

- b.) Fit side dolly through the shoulder bolts and secure assembly with two (2) wing nuts. See [Figure A-22](#).
 - c.) Repeat steps a and b for the other side dolly.

Potential for injury if covers removed and power is left "ON".

- 2.) Disconnect cables on the right side of the rear cover.
- 3.) Remove rear cover.
 - a.) Disengage upper and lower cantrell brackets on both sides of the rear cover.
 - 1.) Using steady but firm pressure, lift each of the lower cantrell brackets from their associated retainers. See [Figure A-9](#).
 - 2.) Disengage the locking mechanism on the upper cantrell brackets by using your thumb to slide the trigger (red lever) back. This will release the locking mechanism and allow the cantrell to be rotated upwards with steady and firm pressure.
 - b.) Disengage the rubber retaining straps on both sides.

CAUTION



4.2 Installation

- 1.) Position cover in back of gantry
- 2.) Attach the rear cover
 - a.) Align the studs on both sides of the rear cover with the receivers located on the gantry frame.
 - b.) Insert the stud on one side into its associated receiver and attach the rubber retaining straps. Then insert the stud on the other side into its associated receiver and attach its rubber retaining straps.

Note: You may find it helpful to lift “up” on the cover to align the stud while attaching the rubber retaining straps.

- 3.) Reattach upper and lower cantrell brackets on both sides.
 - a.) Remove upper cantrell brackets from service position and rotate them into position over their associated retaining pins. Press down firmly on the bracket and snap it into place. The locking mechanism on each upper bracket should lock the bracket securely into place. Do this on both sides.
 - b.) Remove lower cantrell brackets from service position and rotate them into position over their associated retaining pins. Press down firmly on the bracket and snap it into place.

Note: Adjustment of the cantrell brackets can cause misalignment of the top and side covers. The upper and lower cantrell brackets do not re-quire adjust during normal use.

- 4.) Remove dolly, disassemble and store safely away.
- 5.) Re-attach cables to cover.
- 6.) Re-install the mylar (scan) window. Carefully bend the scan window and place it into the channel (groove) provided in the covers.

Figure A-23 Installing the mylar window



Section 5.0

Gantry Cover Alignment Adjustment Guide

5.1 Overview

This section explains the adjustment capability available for the Gantry covers and guidelines for performing adjustments.

5.2 Tools Required


3mm, 5mm, 8mm, Hex wrenches

5.3 Adjustment Guide


5.3.1 Front/Rear Cover adjustment

The front and rear covers are held in a fixed vertical position by the cover locking mechanisms at the vertical middle of each cover on each side. The only adjustment capability is found on each of the cover cantrells (see [Figure A-24](#), one near each corner of the front and rear covers. Adjustments are best done at zero tilt.

Figure A-24 Cantrells

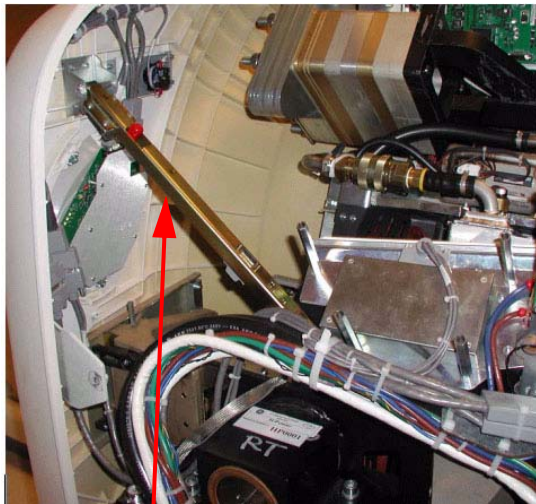

Top Cantrell same
for all top latches

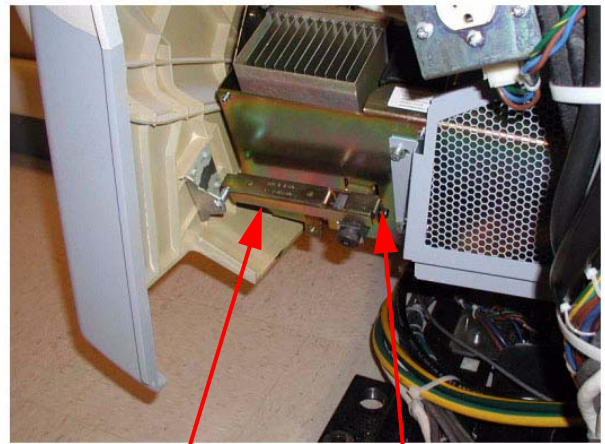

Bottom Cantrell same
for all bottom latches

Adjustment screw on
end of all brackets

- 1.) Remove the Gantry side covers using a 8mm hex wrench to release the bottom latches. Refer to [Section 1.0 on page 113](#), for further details.
- 2.) Disable gantry rotation, 120VAC and HVDC via the service switch panel.
- 3.) Disconnect the fan power via the white molex connector (one for each top cover). Remove the top covers by lifting up and pulling to the side of the gantry.
- 4.) Remove the scan window to avoid creases or kinks during cover movements.
- 5.) Using a 5mm hex wrench, adjust the bolt on the gantry end of the cover cantrells (see [Figure A-24](#)) to lengthen or shorten the cantrell arm that will push or pull the cover corner. The distance

between the inside edges of the front and rear covers will be typically 27" + 1/8" (689mm). Distance should remain the same as measured along the inside edges from top to bottom with both covers hanging level as seen with a level on the inside vertical edge.


5.3.2 Scan Window

Position the scan window in the opening to verify that it seats all the way around with no gaps. The scan window should seat easily if the front and rear covers have been adjusted properly. If not, go back and readjust such that the scan window seats properly. There is no easy solution other than adjusting and looking for fit.

5.3.3 Top Cover Adjustment

The top covers rest in a top bracket that is mounted to the stationary gantry frame and subsequently just rest on top of the front/rear covers. There are also two alignment fingers on the outer edge of the top covers that need to sit in brackets on the front/rear covers.

Figure A-25 Top Cover Bracket

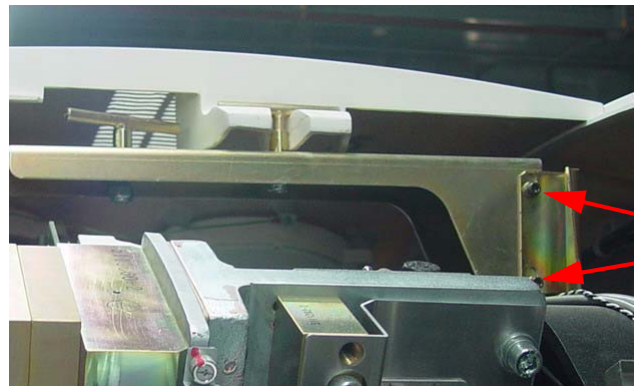

Screws on both sides
allow adjustment up/
down only

- 1.) Set the top covers back on the gantry and look for proper fit. The top covers should rest on the front/rear covers all along the edges. If they do then skip the rest of this section.
- 2.) If the top covers do not rest on the gantry at the top then lift them back off and use a 5mm hex wrench to loosen the bolts that hold the top cover bracket for adjusting vertical alignment to allow the top covers to rest on the front and rear covers.
- 3.) The brackets on the front/rear covers for the fingers on the outer edge of the top covers are adjustable using a 3mm hex wrench and sliding the bracket till it lines up with the top cover finger. (See [Figure A-26](#))

Figure A-26 Top Cover Fingers



- 4.) Verify the top cover tapered edges fit over the lip along the top edge of the front/rear covers. If

not, go back to Section 5.3.1 and readjust the cantrells as necessary.

- 5.) When done adjusting the top covers, reconnect the fan power molex connections.

5.3.4 Side Covers

The side covers just hang from the top covers by two fingers that are not adjustable. When set into the top covers, the side covers should hang straight down with the tapered edges fitting over top of the lip of the front/rear covers.

5.3.5 Finish

- 1.) Re-Enable gantry rotation via the Axial Enable switch on the service switch panel.
- 2.) Replace the side covers and latch using a 8mm hex wrench on the side cover bottom latches.
- 3.) Run an axial scan to ensure proper system operation prior to turning system back over to the customer.

Section 6.0 Gantry Scan Window

6.1 Remove Scan Window

- 1.) Grab the window at the top and pull firmly downward.
- 2.) Continue to pull until the top of the scan window makes contact with the bottom portion of the scan window.
- 3.) Hold the top and bottom portions of the scan window together, grasp both sides of the scan window, move them together and lightly pull upward, until you can free the window from between the front and rear covers.

Figure A-27 Scan Window Removal



6.2 Install Scan Window

- 1.) Install the front and rear covers.
- 2.) Deform the scan window, as shown in [Figure A-28](#), and nest the scan window at the bottom of the opening between the front and rear covers, ([Figure A-29](#)) with the rivets in the 6 o'clock position installation position. Remember the rivets must be in the 12 o'clock position when the mylar window is fully installed.
- 3.) After you complete the initial seating of scan window, let the window slowly unfold, and work both sides of the window into position, starting at the bottom and finishing at the top.
- 4.) Make sure you position the window with the rivets at the 12 o'clock position, and the mylar window slit at either the 3 or 9 o'clock position.

Figure A-28 Install Scan Window


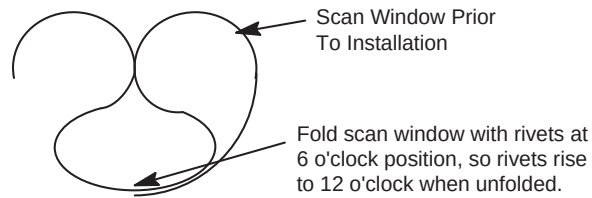


Figure A-29 Scan Window Nested Between Front and Rear Cover

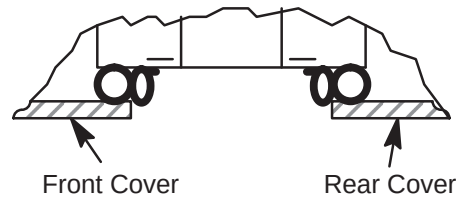

Section 7.0 Gantry Base Covers

Refer to [Figure A-30](#) for the following assembly sequence.

Note: Tighten means torque to 2.3 Nm

- 1.) Position cover (item 5) on gantry base with bracket slots aligned to gantry holes. Center cover left to right and attach with (4) hardware (items 1, 2, 3) as shown and tighten.
- 2.) Assemble two (2) bulkheads (item 14) to two (2) brackets (item 13) using (4) hardware items 3 and 4. Assemble two (2) brackets (item 10) and two (2) brackets (item 13) to gantry base using eight (8) hardware items 1, 2 and 3. Finger tighten hardware with bracket moved outward to end slots. Install side covers (item 6 and 7) on base pushing brackets (item 11 and 13) inward until properly aligned with front cover. Remove side covers, tighten fasteners and replace side covers using one (1) hardware (item 16, 2 and 3) and two (2) item 15 on each cover and tighten.
- 3.) Assemble last bracket (item 11) loosely to gantry base with two (2) hardware (item 1, 2 and 3). Install rear cover (item 8) to base properly aligned to side cover (item 6). Attach rear cover to bracket with hardware (items 1, 2 and 3) tightening all fasteners. Lock latch.
- 4.) Place cover (item 9) on gantry base, aligned to covers (item 7 and 8). Lock both latches.

Figure A-30 Gantry Base Covers


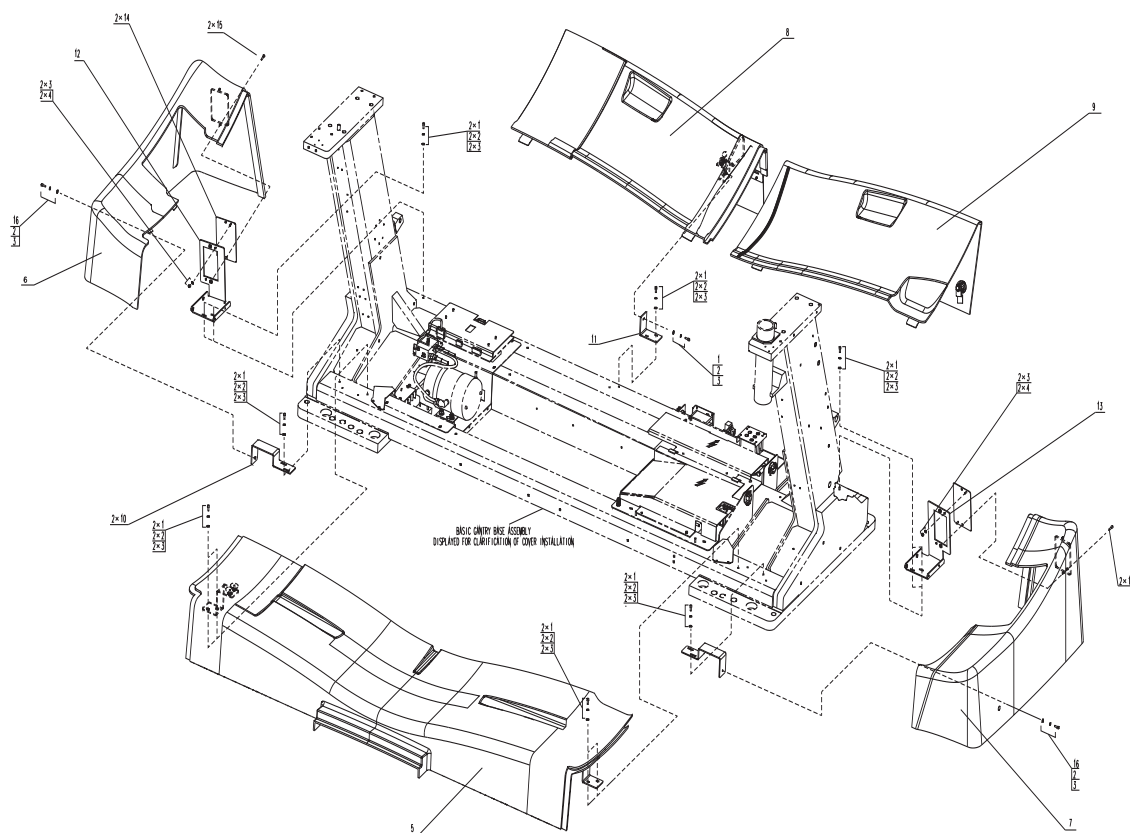


Figure A-32 Installing the Table Flapper

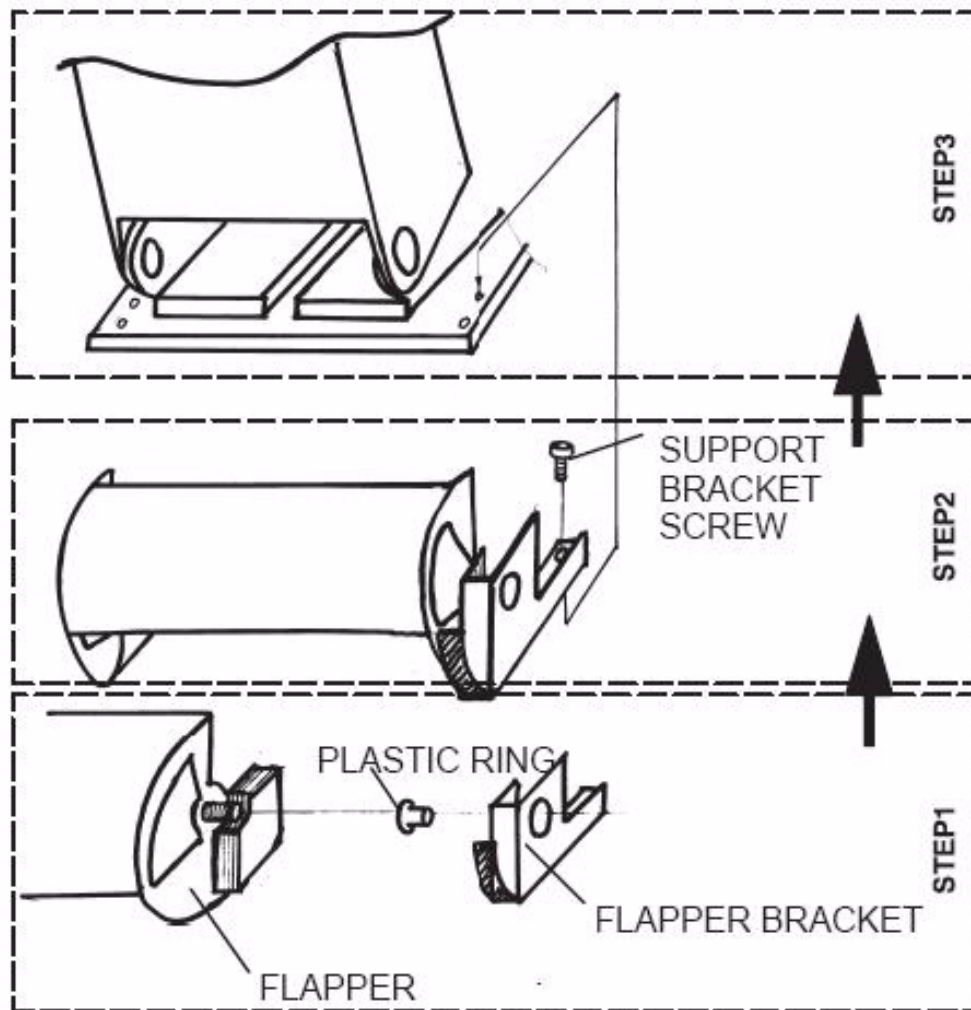

8.3 Skirt Cover installation

- 1.) Attach Skirt Cover to flapper brackets.
- 2.) Mount Skirt Cover to table bracket with two screws.


Figure A-33 Installing Skirt Cover



8.4 Bottom cover installation

- 1.) Install the front bottom cover using four screws. After installation, install the screw caps onto the screws.
- 2.) Install the foot switch pedals using four socket screws. After installation, install the screw caps onto the screws.
- 3.) Install the rear bottom cover so that its tabs are aligned with the frame pins.
- 4.) Pull the cover rearwards, then mount it using two screws.
- 5.) After installation, install the screw caps onto the screws.

Figure A-34 Installing the Bottom Covers

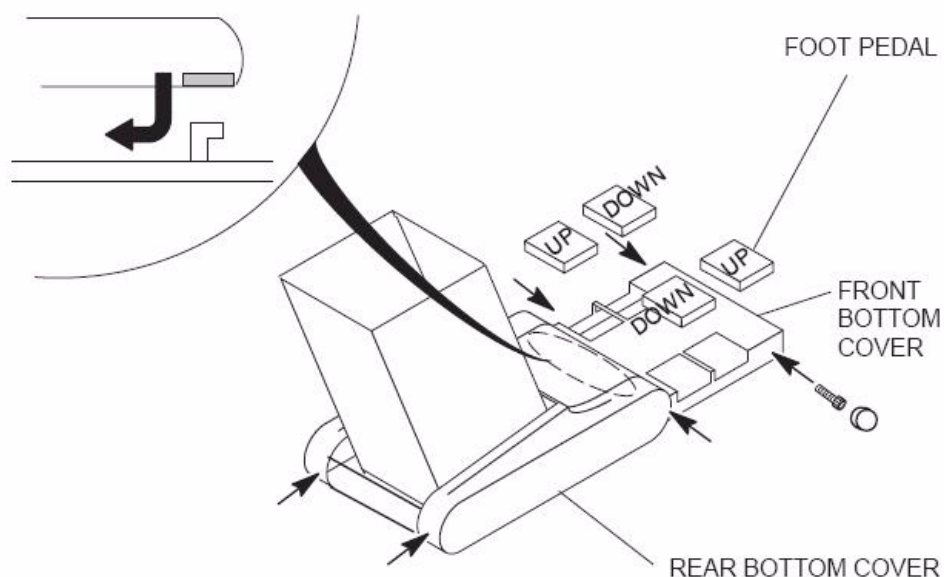

Appendix B

FWS Assembly and Adjustment

Section 1.0

Introduction

Freedom WorkSpace (FWS) is an optional console table with better ergonomic performance. The monitor arm is designed to adjust quickly and easily.

There are two FWS: part 5160735 and part 5168666. They are different in height adjustment. Refer to respective install procedure during installation.

Before assembly go through this section and checklist with the FWS to have an overview.

- Assemble workspace
 - Assemble chain channel and chain with table base (only applicable for part 5160735)
 - Assemble table base and workspace
 - Assemble drawer
- Install monitors
 - Mount pole onto table surface
 - Install monitor with monitor arms
 - Install monitor arms on the pole
 - Route cables
- Adjust the monitors for customer use
- Install Seismic Kit

Figure B-1 Overview of Assembled FWS w/o Peripherals




Section 2.0 Procedure

2.1 Assemble Worksurface (For Part 5160735)


- 1.) Carefully remove the base assembly from the packaging material. With the leg uprights supported, move the legs apart to attach the crossbrace.
- 2.) Secure the crossbrace to the leg uprights with four 1/4 - 20 x 5/16" screws (Figure B-2).
- 3.) Replace the back outer-leg covers to the leg uprights (Figure B-2).
 - a.) Remove the adhesive pad paper on the legs.
 - b.) Position the top of each cover up into the upright grommet.
 - c.) Press the bottom down to secure with the adhesive pad.

Figure B-2

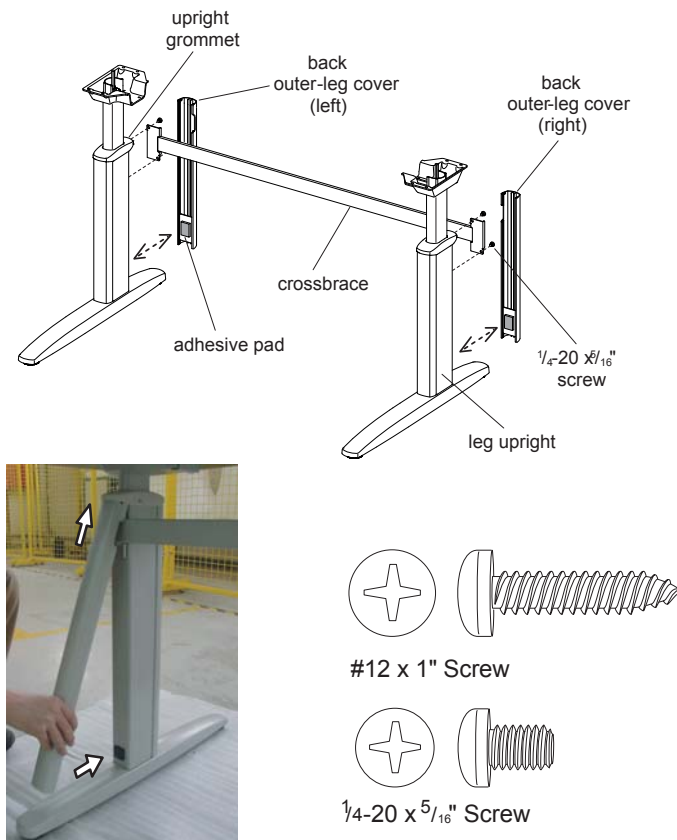


- 4.) Slide the chain channel into the two openings in the worksurface supports ([Figure B-3](#)).

Figure B-3

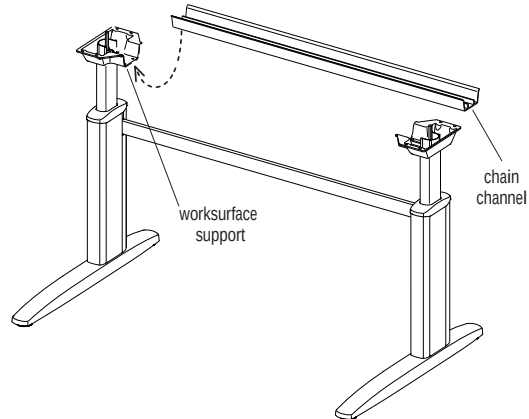

- 5.) Insert hand crank on one hex shaft, adjust the legs height and ensure the chain channel is horizontal.
- 6.) Install chain ([Figure B-4](#)).
 - a.) Pull out the cotter pin.
 - b.) Slide the sprocket off over the hex shaft.
 - c.) Engage the sprocket into the chain and reinstall the sprocket.
 - d.) Insert the cotter pin and spread the ends to secure.
 - e.) Repeat the steps on the other side.
 - f.) Insert hand crank on one hex shaft and rotate to ensure the worksurface is adjustable.

Figure B-4

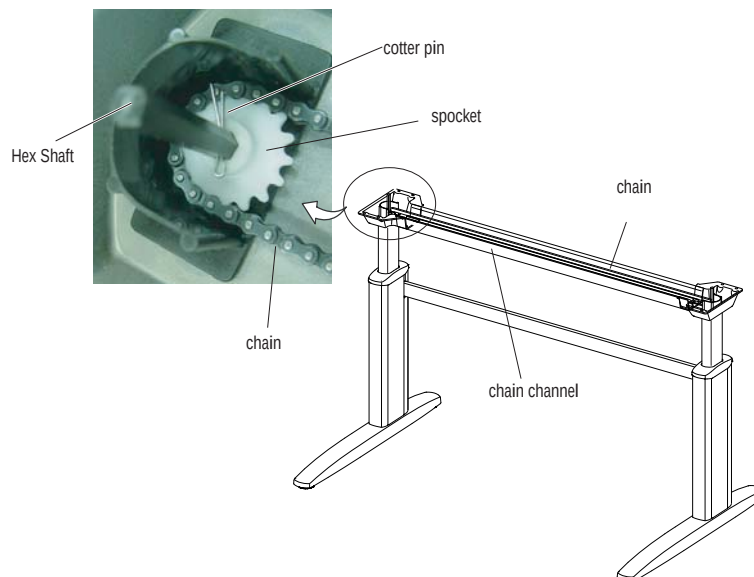


- 7.) Place the worksurface upside down on a clean, soft surface to prevent scratching. Invert the table base onto the worksurface, taking care to align the holes with the pre-drilled holes in the worksurface. Secure the table base to the worksurface with eight #12 x 1" screws ([Figure B-5](#)).

Figure B-5




- 8.) Install two Chain Channel Brackets using four metal screws and two plastic push screws. (Figure B-5)
- 9.) Find the 8 drawer screws inside the drawer. Install the drawer. (Figure B-6)

Figure B-6

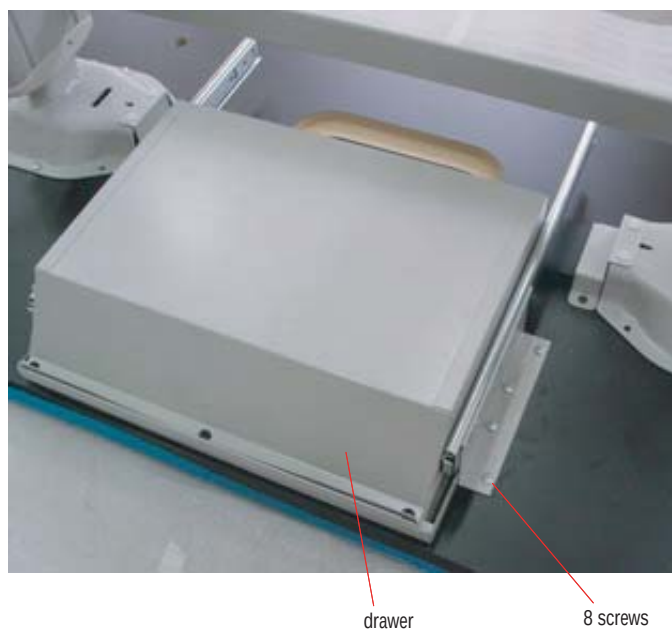

- 10.) Carefully turn the table assembly to the upright position. Insert crank grommets into holes in the worksurface as illustrated. Insert crank into one grommet and place a filler cap into the other side. Adjust worksurface to desired height with hand crank (Figure B-7).

Figure B-7

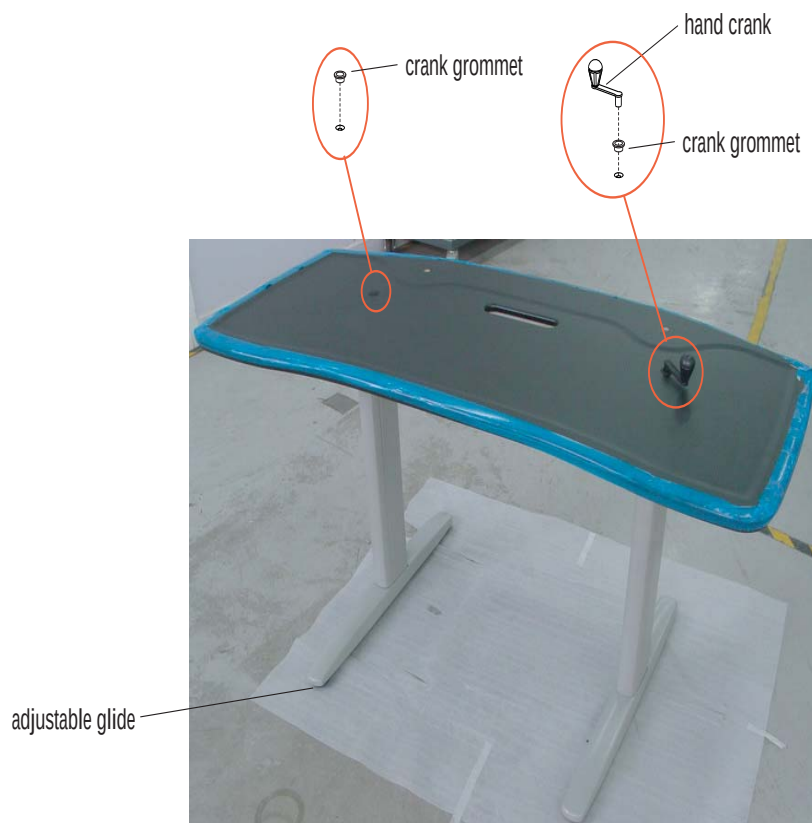

- 11.) The four adjustable glides may be raised or lowered to accommodate varying floor surfaces. Lift the leg off the floor and turn the glide clockwise or counterclockwise to adjust (Figure B-7).

2.2 Assemble Worksurface (For Part 5168666)

- 1.) Carefully remove the base assembly from the packaging material.
- 2.) Attach the crossbrace with the two inner-leg cover and secure the crossbrace to the down leg with three M6 x 20 screws (Figure B-9).

Figure B-8

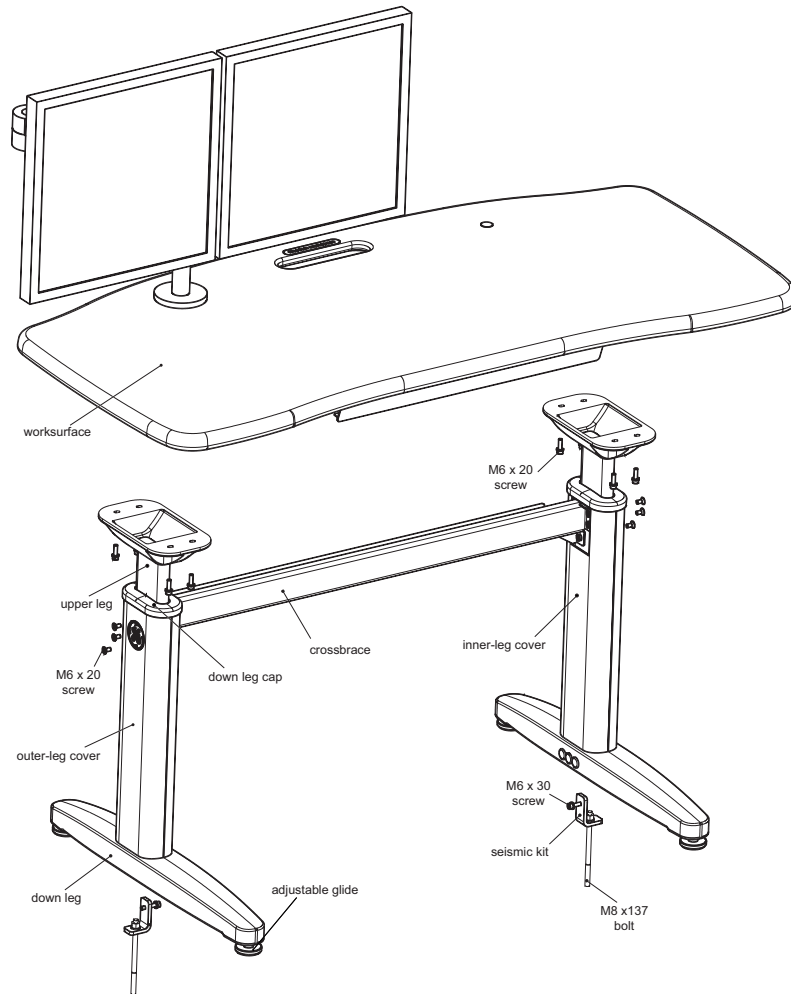


Figure B-9



- 3.) Press the inner-leg cover to secure with the adhesive pad ([Figure B-10](#)).

Figure B-10



- 4.) Install upper legs to down legs.
 - a.) Install down leg cap to the down legs ([Figure B-11](#)).
 - b.) Insert upper legs into down legs, adjust to desired height by attaching the upper legs' screws to appropriate down legs' holes. [Figure B-12](#) shows an example of attaching upper legs at the fifth and seventh holes from the bottom upward. In this case, the height of worksurface would be 785mm.
 - c.) Secure the upper legs with four M6 x 20 screws each leg.

Figure B-11



Figure B-12



- 5.) Insert outer-leg cover and inner-leg cover into down leg slot. Then push down the upper cover to fix them.([Figure B-13](#)).

Figure B-13



- 6.) Place the worksurface upside down on a clean, soft surface to prevent scratching. Invert the table base onto the worksurface, taking care to align the holes with the pre-drilled holes in the worksurface. Secure the table base to the worksurface with eight M6 x 20 screws ([Figure B-14](#)).

Figure B-14



- 7.) Carefully turn the table assembly to the upright position. Insert hole cover into one of the two holes in the worksurface, the other hole would be used for monitor installation.
- 8.) The four adjustable glides may be raised or lowered to accommodate varying floor surfaces. Lift the leg off the floor and turn the glide clockwise or counterclockwise to adjust.

2.3 Install Monitors

Note: For more information, refer to the materials shipped with the FWS.

Figure B-15 (Part 5160735)

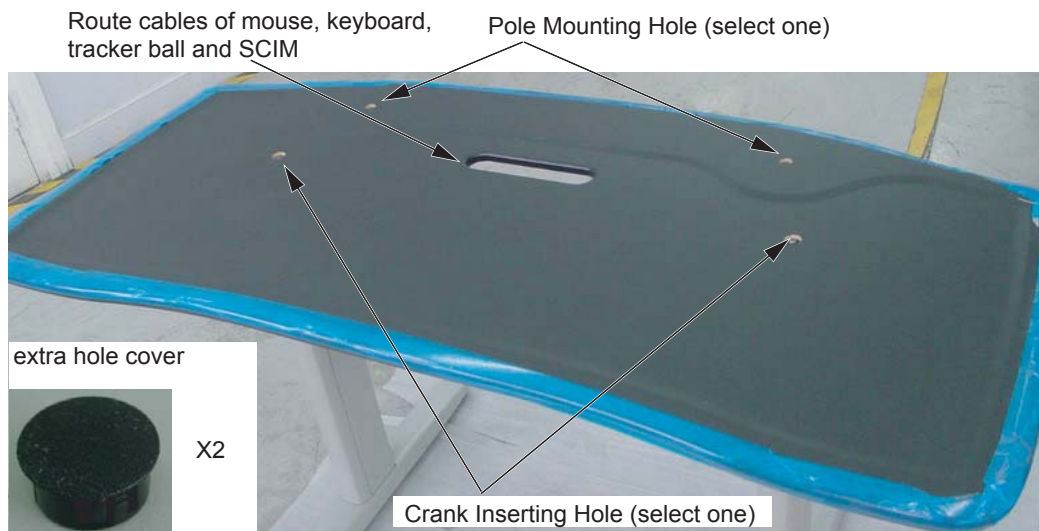



Figure B-16 (Part 5168666)



- 1.) Mount the pole on the worksurface. (See [Figure B-15](#), [Figure B-16](#) and [Figure B-17](#))
 - a.) Adhesive the pad under the grommet mount.
 - b.) Use large bolt to go through parts.
 - c.) Fix the clamp on the desk with 4 screws.

d.) Put the screw plate cover through the pole. Check the pin position with the screw.

Figure B-17

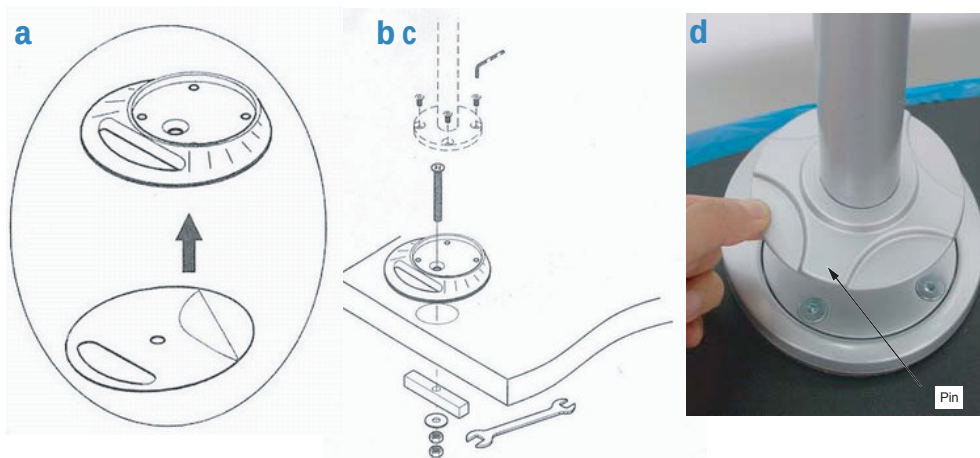


- 2.) Change configuration of the monitor arms for site.
 - a.) Three configurations for customer to select, see pictures below:

Figure B-18



- b.) If B or C is selected, disengage the not-to-used arm by unscrewing the junction bolts and reconfigure the arms.

Figure B-19

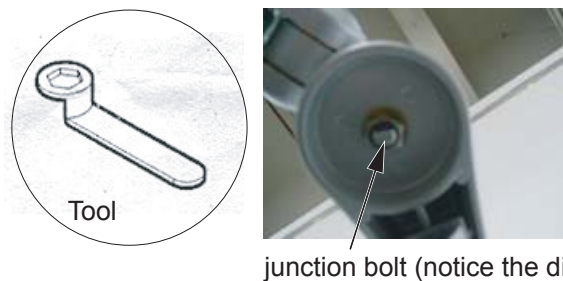


- 3.) Secure cable covers onto monitor arms. Cable covers should open backward the customer.

Figure B-20




- 4.) Mount the monitors to the arm.

Figure B-21

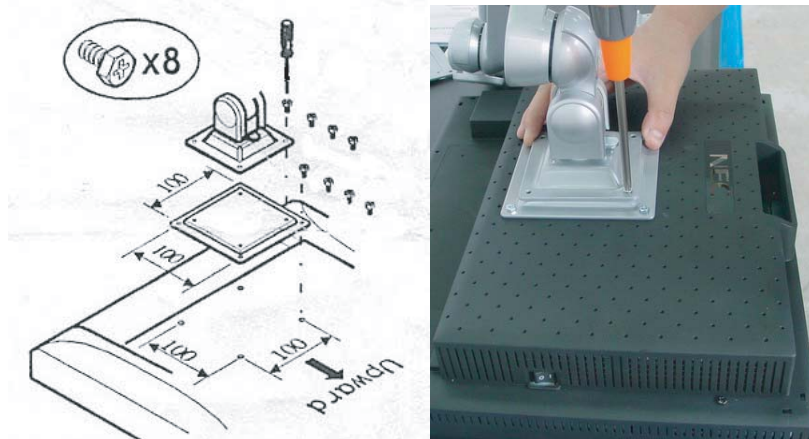


- 5.) Mount arms with pole. (Figure B-22)
- Loosen two screws slightly. Put the arm into the pole and adjust it to the suitable position. Middle joint mount (two large allen screws) of the monitor arms should face toward the customer.
 - Tighten up the pole with 4mm and 3mm allen keys.
 - Put in the plug on top of the pole.

Figure B-22

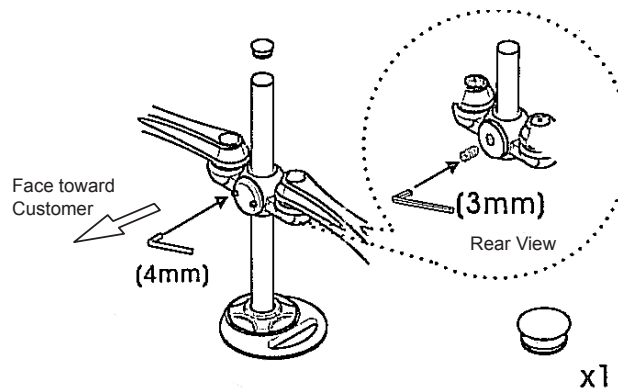

- 6.) Put extra hole covers onto the holes not used. (Figure B-15)
- 7.) Route cables.
- Extended Cables Kit (5160577) for FWS is used to connect monitors, keyboard, mouse, SCIM and trackball. Connection is same with short cables, see [Section 4.0 on page 79](#), Console Connections.
 - Thread monitor cables through cable covers.
 - Use cable tie to wrap the cables together or wrap the cables with pole.

Figure B-23

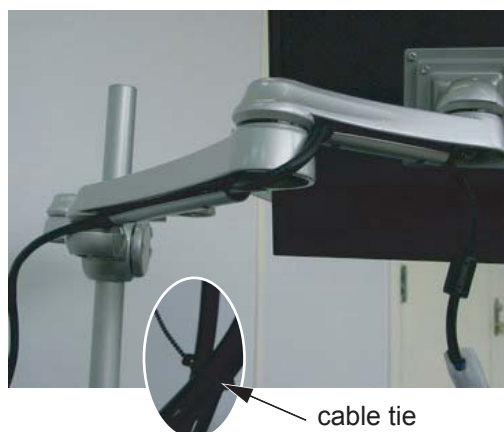

- 8.) Route cables (for FWS Part 5168666), thread cables through worksurface and use cable tie to wrap them together to the back slot of crossbrace ([Figure B-24](#)).

Figure B-24




2.4 Monitor Adjustment

- 1.) Position the arms for ergonomic viewing. For optimum ergonomic viewing set top of monitor screen 1" below eye height. With monitor and arm properly supported, loosen screw and adjust up or down to desired height. Retighten to secure arm to pole.
- 2.) Adjust monitor tilt: Tilt monitor back and forth through entire pivot range of motion. If it does not stay in place or movement in one direction is stiff, pivot needs adjustment - see below (pivot cover is removed). Adjust until monitor stays in place and tilting forces, up or down, are equal.
Increase tilt lifting force: Turn screw clockwise.
Decrease tilt lifting force: Turn screw counterclockwise.

Figure B-25

side view

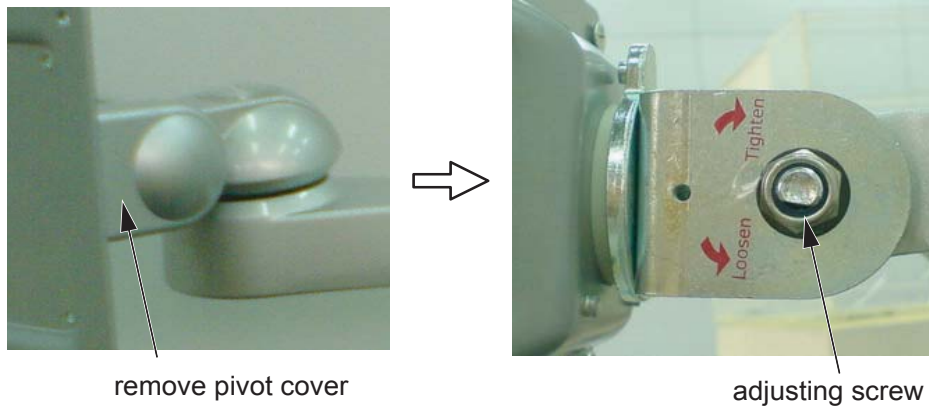

2.5 Install Seismic Kit (For Part 5160735)

If site specifications require seismic mounting, follow below steps:

- 1.) Insert the seismic bracket (x4) into the gap above the adjustable glide, at the bottom of leg.
- 2.) Use 1/4" bolts to mount the brackets to the floor. The seismic brackets are shipped with FWS table.

Figure B-26



2.6 Install Seismic Kit (For Part 5168666)

If site specifications require seismic mounting, follow below steps:

- 1.) Use M8 x 137 bolts to mount the brackets to the floor at a distance 916mm (Figure B-27).
- 2.) Mount table base to brackets with two M6 x 30 screws (Figure B-28).

Figure B-27


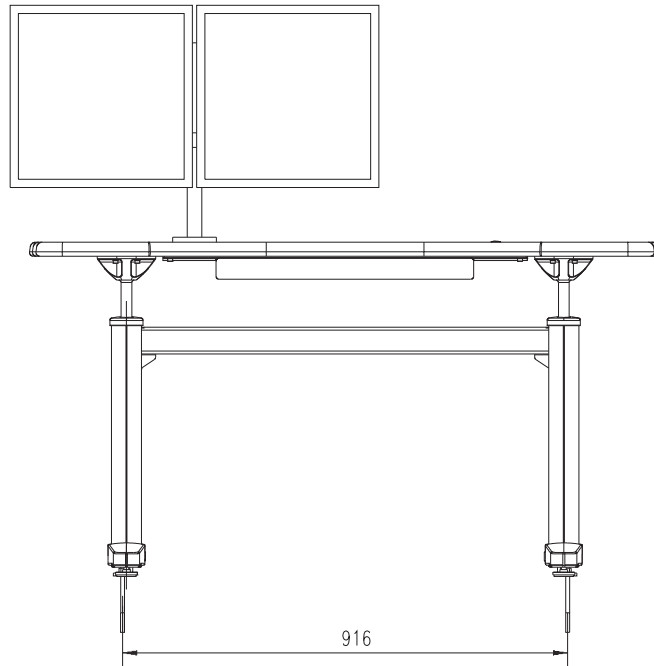


Figure B-28

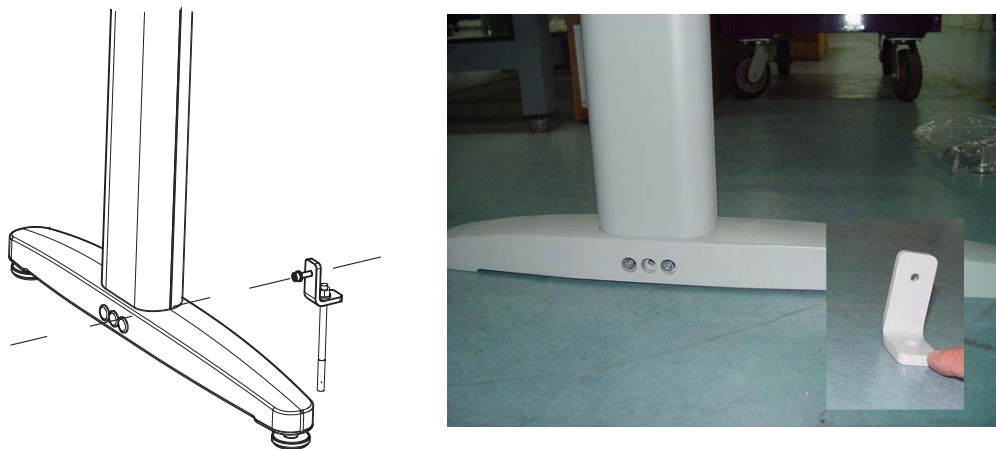

Appendix C

Pictorial Representation of Required Tools

Use the following guide as a reference, if you are unsure of a tool listed in [Section 2.5, on page 26](#).

Table C-1 Required Tools

TOOL NAME	PICTURE	EXAMPLE PART NUMBER*
Adapter		Sears Industrial: $\frac{3}{8}$ " to $\frac{1}{2}$ " (9-4258)
Ball-Peen Hammer		Sears Industrial: 1lb/2lb (9-38465)
Canned Air		Miller Stephenson: Aero Duster (MS-222N)
Clamp on Amp Meter		Sears Industrial: 9-WTAD105
Combination Wrench Set		Sears Industrial: U.S. Standard & Metric (9-44048)
Cordless Screwdriver		Sears Industrial: 9-MU65401
Deep Well Socket		Sears Industrial: $\frac{3}{4}$ " X $\frac{3}{8}$ " (included with 9-34496)
Dental Pick		
Diagonal Cutting Pliers		Sears Industrial: Small (9-45077)

* Part Numbers given for reference only. GEMs does not endorse any tool brand name.

Table C-1 Required Tools (Continued)

TOOL NAME	PICTURE	EXAMPLE PART NUMBER*
Drill		Sears Industrial: $\frac{3}{8}$ " or $\frac{1}{2}$ " (9-27859)
Drill Adapter		Sears Industrial: 3" X $\frac{3}{8}$ " (9-APSZ24)
Drill Bit Set		Sears Industrial: U.S. Standard (9-66084)
DVM		Sears Industrial: 9-82028 Sears Industrial: 9-FL873
Extension for Ratchet Wrench		Sears Industrial: 3" X $\frac{1}{2}$ " (9-44133)
Gloves		Sears Industrial: Large (9-40502)
Hammer Drill		Sears Industrial: $\frac{1}{2}$ " (9-27205)
Hex Bit Set		Sears Industrial: $\frac{1}{4}$ " (9-SK45508)
Hex Key (Allen Wrench) Set		Sears Industrial: U.S. Standard (9-46284)

* Part Numbers given for reference only. GEMs does not endorse any tool brand name.

Table C-1 Required Tools (Continued)

TOOL NAME	PICTURE	EXAMPLE PART NUMBER*
Level		Sears Industrial: 4' (9-39856)
Masonry Bit		
Open-End Wrench (Thin or Standard Tappet)		Snap-on: 10mm (SRSM10) & 21mm (LTAM2124)
Ratchet Wrench		Sears Industrial: $\frac{3}{8}$ " (9-43175)
Reciprocating Saw with Blades		Sears Industrial: 9-MU650921
Safety Glasses		Sears Industrial: 9-18650
Safety Shoes		
Screwdriver Set		Sears Industrial: Phillips & Straight (9-41505)
Socket Set		Sears Industrial: Standard $\frac{3}{8}$ " (9-34496)
Sockets		Sears Industrial: 1 $\frac{1}{8}$ " X 1 $\frac{1}{2}$ " (9-47516)

* Part Numbers given for reference only. GEMs does not endorse any tool brand name.

Table C-1 Required Tools (Continued)

TOOL NAME	PICTURE	EXAMPLE PART NUMBER*
Step Ladder		Sears Industrial: 6' (9-WN6006)
Tongue & Groove Pliers		Sears Industrial: Large (9-CL440)
Torpedo Level		Sears Industrial: 9" (9-39829)
Torque Wrench		Sears Industrial: $\frac{3}{8}$ " (9-WR3470)
Universal Joint		Sears Industrial: $\frac{3}{8}$ " (9-4435)
Vacuum Cleaner		Sears Industrial: 8 Gal (9-17780)
* Part Numbers given for reference only. GEMs does not endorse any tool brand name.		

AppendixD

Option Cables

Run all option cables that connect to the Table, Gantry or Console.
For example the remote monitor (through the splitter) and SmartScore.
For more information, refer to “LightSpeed Linux CRT & LCD Monitor Option”

Figure D-1 Remote Monitor Splitter



Figure D-2 SmartScore



AppendixE

Regulatory Clearance Quick Reference Guide

Section 1.0

Regulatory Code Description

Egress: 29 CFR 1910 Subpart E (OSHA) and NFPA 101 (Life Safety Code) define the minimum requirements for means of egress. The requirement most applicable to equipment installation and room layout is minimum width of exit access. Under OSHA 1910.37(f)(6), the minimum width of exit access shall in no case be less than 28 in. from any potentially occupied point in the room.

Under NFPA 101 (2006 edition) 7.3.4.1, the minimum width of any means of egress is 36 in. However, NFPA allows this to be reduced to 28 in. around furniture or equipment, provided that a 36 in. clearance would otherwise be available without moving permanent walls.

Electrical Clearance: 29 CFR 1910 Subpart S (OSHA) and NFPA 70E (Standard for Electrical Safety in the Workplace) define minimum clearance requirements for the workspace around electrical equipment. Under both OSHA 1910.303(g)(1) and NFPA 70E (2004 edition) 400.15, a minimum clear space of 36 in. depth (with minimum 30 in. width and 78 in. height) must be provided in front of electrical equipment with parts operating at 600 volts or below and likely to require examination, adjustment, servicing, or maintenance while energized.

This safety clearance requirement applies to all GEHC equipment. Although 36 in. is the minimum clearance for most installations, the standards require an increased minimum clearance distance where parts operate above 150 volts (but still below 600 volts) under the following circumstances:

- If the wall or surface directly facing the electrical equipment is grounded (e.g. brick, concrete, or tile) or includes grounded protrusions (such as medical gas ports, metal door or window frames, water sources and metallic sink structures, metallic cabinetry, electrical disconnects or emergency off panels, air conditioners or vents), then a 42 in. clearance depth is required.
- If the possibility exists of exposed and unguarded live parts on both sides of the workspace (for example if a power distribution unit were positioned on the wall directly facing the GEHC equipment), then a 48 in. clearance depth is required.

Section 2.0

Terms and Definitions

Egress: The path of exit from within any room. U.S. regulations require a minimum of 28 in. (711.2 mm) of continuous and unobstructed space, including trip hazards along the path of exit.

Workspace: The dimensional box required for safe inspection or service of energized equipment. It consists of depth, width, and height. The depth dimension is measured perpendicular to the direction of access. The U.S. regulation minimum is 36 in. (914.4 mm), but additional conditions can increase the minimum dimension requirement. GEHC defines this as the envelope of the component superstructure with the external covers in place.

Service Access Width: The width of the workspace in front of the equipment. A minimum of 30 in. (762 mm), or the width of the equipment, whichever is greater.

Head Clearance: The height dimension of the workspace. The height of the workspace measured from the floor at the front edge of the equipment to the ceiling or overhead obstruction(s). 78 in. (1981.2 mm), or the height of the equipment, whichever is greater.

Grounded Wall: Any wall that can be electrically conductive to earth ground. Masonry, concrete, and tile are considered conductive. Additional commonly found aspects of a wall should also be considered grounded. This is not an all-inclusive list:

- Medical gas ports and plates
- Metal doors and window frames
- Water sources and metallic sink structures
- Metallic wall-mounted cabinetry
- A1 main disconnect panel
- Equipment Emergency Off panels
- Industrial equipment (such as air conditioners and vents)
- Expansion joints
- Surface raceway
- Exposed wall conduits
- Floor outlets boxes

The following are not considered as grounded elements of a common wall:

- Standard wall outlet
- Light switches
- Telephones
- Communication wall jacks
- Ceiling tile grids

Section 3.0

Regulatory Minimum Working Clearance by Major Subsystem

Requirements apply to equipment operating at 600V or less, where examination, adjustment, servicing, or maintenance is likely to be performed while live parts are exposed.

Direction of Service Access is defined as perpendicular to the surface of the equipment being serviced.

- Required regulatory clearance distances must be maintained and may not be used for storage. This includes normal system operation as well as service inspection or maintenance..

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access: front of console	914 mm (36 in.)	There are no exposed live part hazards with the cover in place. If the console is placed under a counter, the front edge of the console must be even with the vertical edge of the console workspace. Note: This component is typically serviced from the front with access to the rear.
Service access width: Front of console	762 mm (30 in.)	This is the width of the workspace in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.
Head clearance	1981.2 mm (78 in.)	This is the height of the workspace measured from the floor at the front edge of the equipment to the ceiling or overhead obstruction(s). A minimum of 1981.2 mm (78 in.) or the height of the equipment, whichever is greater, is required.

Table 3-8 Console Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (Front of NGPDU)	914.4 mm (36 in.)*	There are no exposed live part hazards with the cover in place. This component is typically serviced from the front with access to the rear. *If exposed live parts of 151 - 600 volts are present, 1219 mm (48 in.) is required on both sides of the workspace with the operator between. *If the opposite wall is grounded and exposed live parts of 151 - 600 volts are present, 1067 mm (42 in.) is required.

Table 3-9 NGPDU Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Service Access Width (Left-Right of workspace)	762 mm (30 in.)	This is the width of the working space in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.
Head Clearance	1981 mm (78 in.)	This is the height of the workspace measured from the floor at the front edge of the equipment to the ceiling or overhead obstruction(s). A minimum of 1981 mm (78 in.) or the height of the equipment, whichever is greater, is required.

Table 3-9 NGPDU Subsystem

- For the gantry and table, distances are measured from the enclosure, not the finish covers.

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (All Sides)	914 mm (36 in.)	If exposed live parts of 151 - 600 volts are present, 1219 mm (48 in.) on both sides of workspace with the operator between is required. If the opposite wall is grounded and exposed live parts of 151 - 600 volts are present, 1067 mm (42 in.) is required.
Service Access Width (Left-Right of workspace)	762 mm (30 in.)	This is the width of the working space in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.

Table 3-10 Gantry Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (Table Head or Foot)	914 mm (36 in.)	There are no exposed live parts hazards with the cover in place. This component is typically serviced from all four sides. this is the width of the workspace on each side of the equipment. A minimum of 914.4 mm (36 in.), or the width of the equipment, whichever is greater, is required.
Direction of Service Access (Table Sides)	914 mm (36 in.)*	*This distance can be reduced to 711 mm (28 in.) provided a written and signed approval is obtained by the local team from the local AHJ (Authority Having Jurisdiction). The signed document must be on file with GE.

Table 3-11 Table Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (Table Foot)	711 mm (28 in.)	For the front gantry cover removal, a minimum of 457 mm (18 in.) is allowed only if an unobstructed egress space of 711 mm (28 in.) is maintained around the equipment for room exit. This also means no trip hazards exist along the path of egress.
Service Access Width (Left-Right of workspace)	762 mm (30 in.)	This is the width of the working space in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.

Table 3-11 Table Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (Front of UPS)	914.4 mm (36 in.)*	There are no exposed live part hazards with the cover in place. This component is typically serviced from the front with access to the rear. *If exposed live parts of 151 - 600 volts are present, 1219 mm (48 in.) is required on both sides of the workspace with the operator between. *If the opposite wall is grounded and exposed live parts of 151 - 600 volts are present, 1067 mm (42 in.) is required.
Service Access Width (Right side and length of UPS)	762 mm (30 in.)	This is the width of the working space in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.
Head Clearance	1981 mm (78 in.)	This is the height of the workspace measured from the floor at the front edge of the equipment to the ceiling or overhead obstruction(s). A minimum of 1981 mm (78 in.) or the height of the equipment, whichever is greater, is required.

Table 3-12 UPS Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (Front of A1 Disconnect)	914.4 mm (36 in.)*	There are no exposed live part hazards with the cover in place. This component is typically serviced from the front with access to the rear. *If exposed live parts of 151 - 600 volts are present, 1219 mm (48 in.) is required on both sides of the workspace with the operator between. *If the opposite wall is grounded and exposed live parts of 151 - 600 volts are present, 1067 mm (42 in.) is required.
Service Access Width (Right side and length of A1 Disconnect)	762 mm (30 in.)	This is the width of the working space in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.
Head Clearance	1981 mm (78 in.)	This is the height of the workspace measured from the floor at the front edge of the equipment to the ceiling or overhead obstruction(s). A minimum of 1981 mm (78 in.) or the height of the equipment, whichever is greater, is required.

Table 3-13 A1 Disconnect Subsystem

Section 4.0

Additional Regulatory Clearance Information

The CT Gantry Left Side Limited Access Initiative provides the capability to reduce the minimum room size for CT Systems while still meeting all installation requirements and specifications. This adds left side flexibility, allowing the CT system to be sited in rooms with widths 22 inches (559 mm) smaller than the current minimum room width.

In the minimum room configuration, left-side access and egress may be restricted for sites with less than 28 inches (711 mm) of left side clearance. Refer to your site's installation print for your room's detail.

4.1 Regulatory Caution

Site prints are required for all system installations including relocation and moves. CT room layout as shown on your site print shall meet all regulatory requirements as described in the installation manual. Additional room components such as cabinets reduce room size. Equipment not shown on the site print may void the caution statement, making the room non-compliant. Actual site measurements before installation will be taken to determine room size and compliance.

For minimum room size sites, where the Gantry distance between the wall and side cover is between 14 inches (356 mm) and 27 inches (686 mm) the left side cannot be considered or used for egress. In this configuration, egress around the left side of the Gantry will be restricted. Unobstructed egress route shall be provided around the Table foot end or behind the back of the Gantry. If mill-work is required or will be installed later, consider room placement locations that allows for regulatory compliance and clearances.

4.2 Operational Caution

In this minimum room layout (14 in. to 27 in. or 356 mm to 686 mm) the customer should consider workflow, customer access for patient care, and critical-care operations space requirements. Additionally, there may be limited equipment access on the Gantry left side when loading patients or when positioning patient equipment in the room between the Gantry and the wall. Detailed customer installation tasks are detailed in the product PIM Chapters 1-4.

4.3 Recommended Room Size

This room configuration offers the most flexibility for future upgrades. It has sufficient workspace and space to add millwork, while meeting all regulatory requirements. This room would be compatible with most two-step future installations.

4.4 Typical Room Size

This room configuration allows for some future upgrades. It has sufficient workspace but limited space to add millwork and meet all regulatory requirements. This room may be compatible with some two-step future installations.

4.5 Minimum Room Size

This room configuration allows for no future upgrades. It has limited workspace and no in-room millwork, but meets all regulatory requirements. This room is not compatible with two-step future installations.

Section 5.0 System Specifications (BrightSpeed Select)

Recommended Room Size	Typical Room Size	Minimum Room Size ³
3962 mm x 6299 mm (13 ft. x 20 ft. 8 in.) ¹	3912 mm x 5868 mm (12 ft. 10 in. x 19 ft. 3 in.) ¹	3310 mm x 5918 mm (10 ft. 10 in. x 19 ft. 5 in.) ²

1 Same Regulatory requirements apply.

2 Same Regulatory requirements apply, with the addition of no energized left side service. No energized left side service is standard on all limited access compatible systems.

3 Note: Use 2040mm gantry width for Minimum Room sitting. There is no electrical parts outside of 2040mm.


5.1 Rooms w/ Less Than 711mm (28 in) Egress Clearance at Table Foot End

Egress requires a clear unobstructed route out of the room, either around the back of the gantry or around the foot end of the table. If your egress route is not around the foot end of the table, maintain 457mm (18 in.) of clearance between the foot end of the table and the wall. Additionally, there needs to be a continuous width of 3200mm (126 in.) - 1600mm (63 in.) on each side of the table center line - to any obstruction so that the front cover can be removed. This may limit placement of millwork within the room. Refer to the Pre-Installation Manual for more details on service clearances.



CAUTION


Rooms smaller than 10 ft 10 in. x 19 ft 5 in. (3310mm x 5918mm) may require construction to meet the minimum requirements. The design center or your project manager of installation may have additional recommendations for your room size.

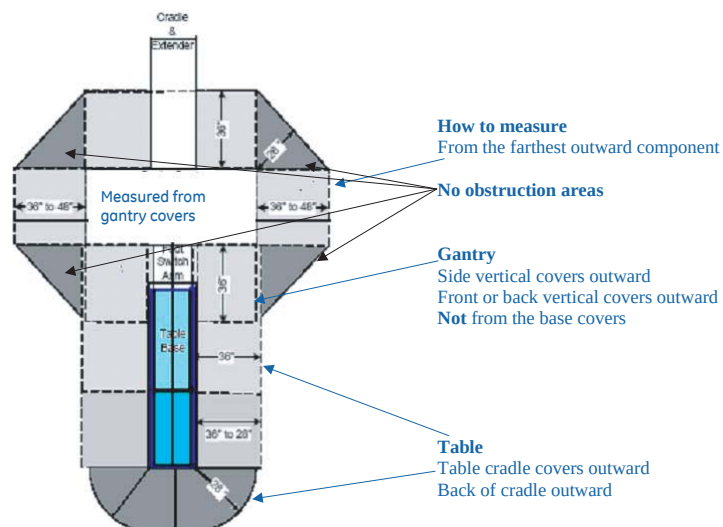


Section 7.0

Minimum Room Size & Requirement Layouts

Room A - Less than 711 mm (28 in.) but greater than 356mm (14 in.) measured from the covers to the left sidewall. In this configuration service, egress and workspace are compromised around the gantry.

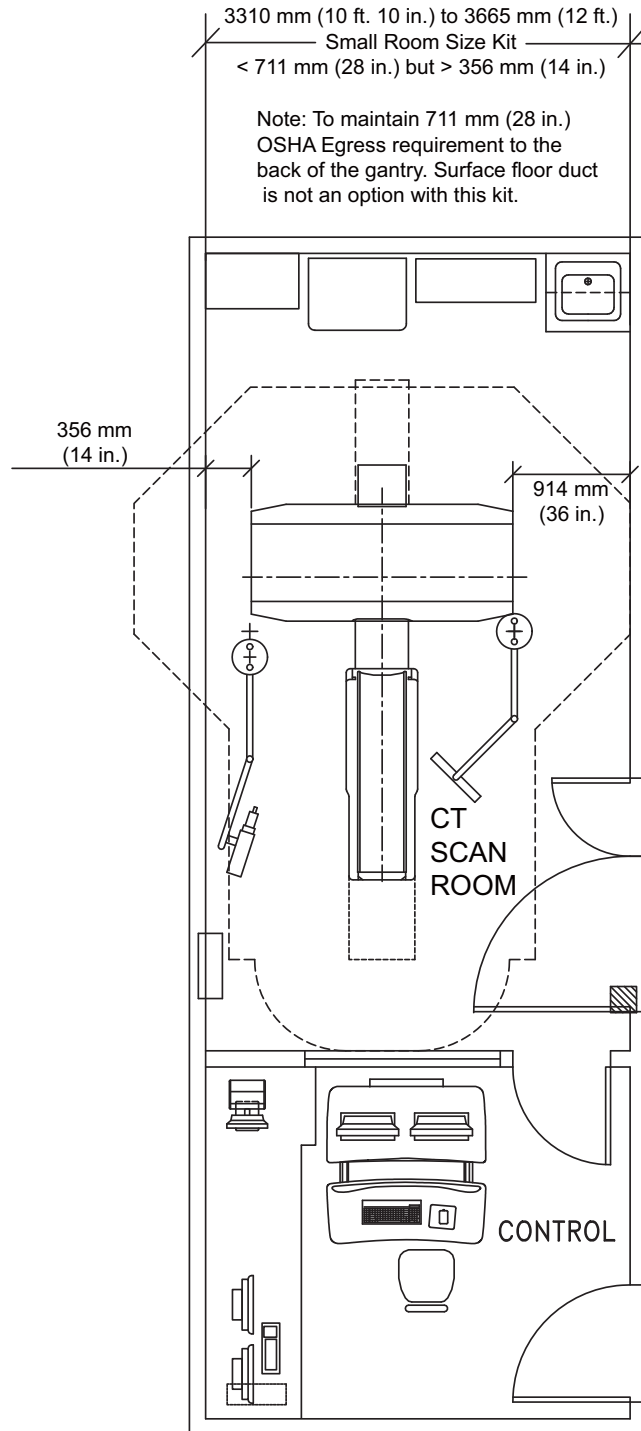

Room B - Less than 914 mm (36 in.) but greater than 711 mm (28 in.) measured from the covers to the left sidewall. In this configuration service, egress and workspace are acceptable around the gantry.

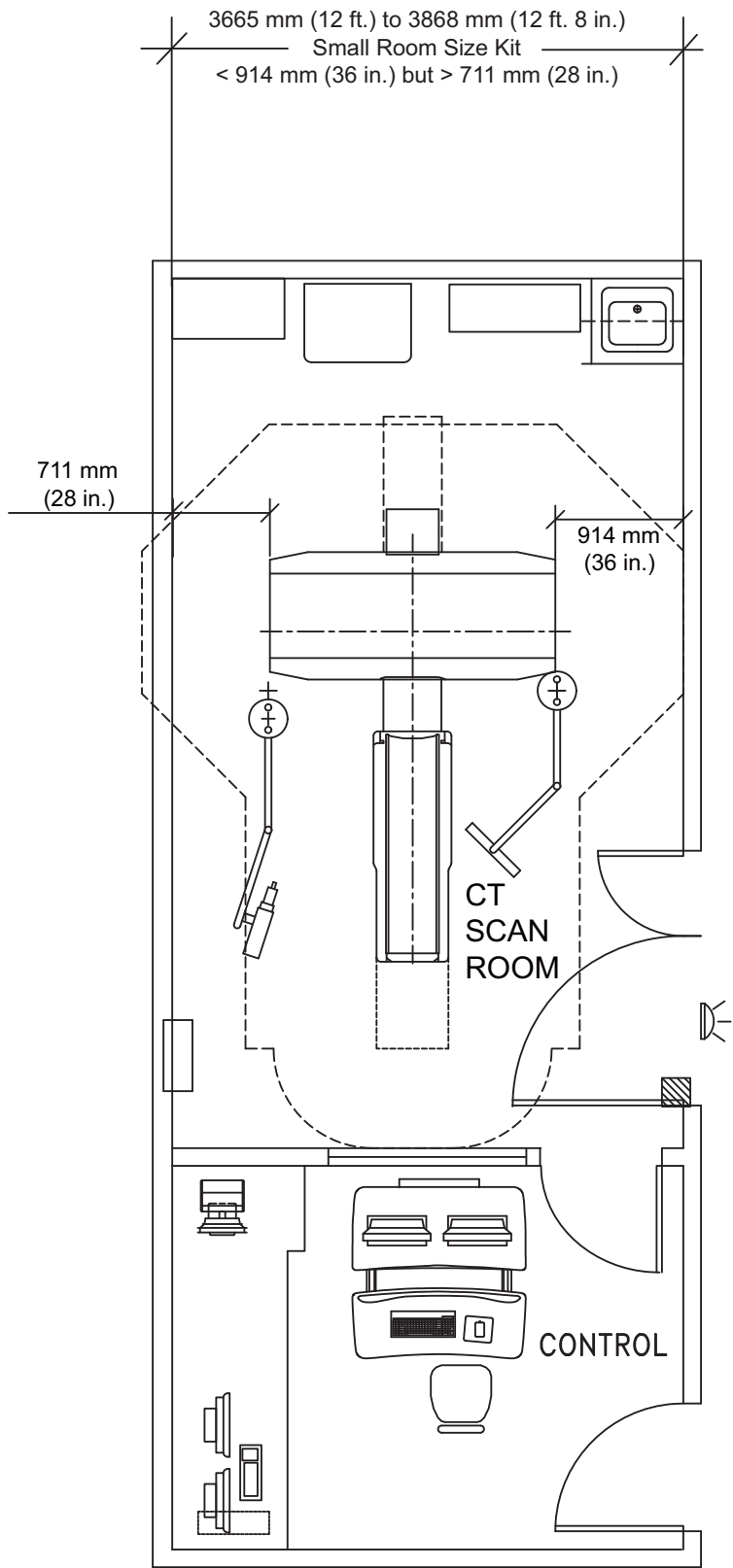

7.1 Typical and Recommended Room Size & Requirement Layouts

Figure E-1, Figure E-2 and Figure E-3 show Typical and Recommended room layouts. You need to be aware of locations for me gas, surface ductwork, or other items that make a grounded wall.

Note: Your room layout may meet the Typical or Recommended room requirements but look different than Figure E-1, Figure E-2 and Figure E-3. Contact your sales person to have a detail room layout completed for your site.

Figure E-1 Typical Room Layout for BrightSpeed Select System

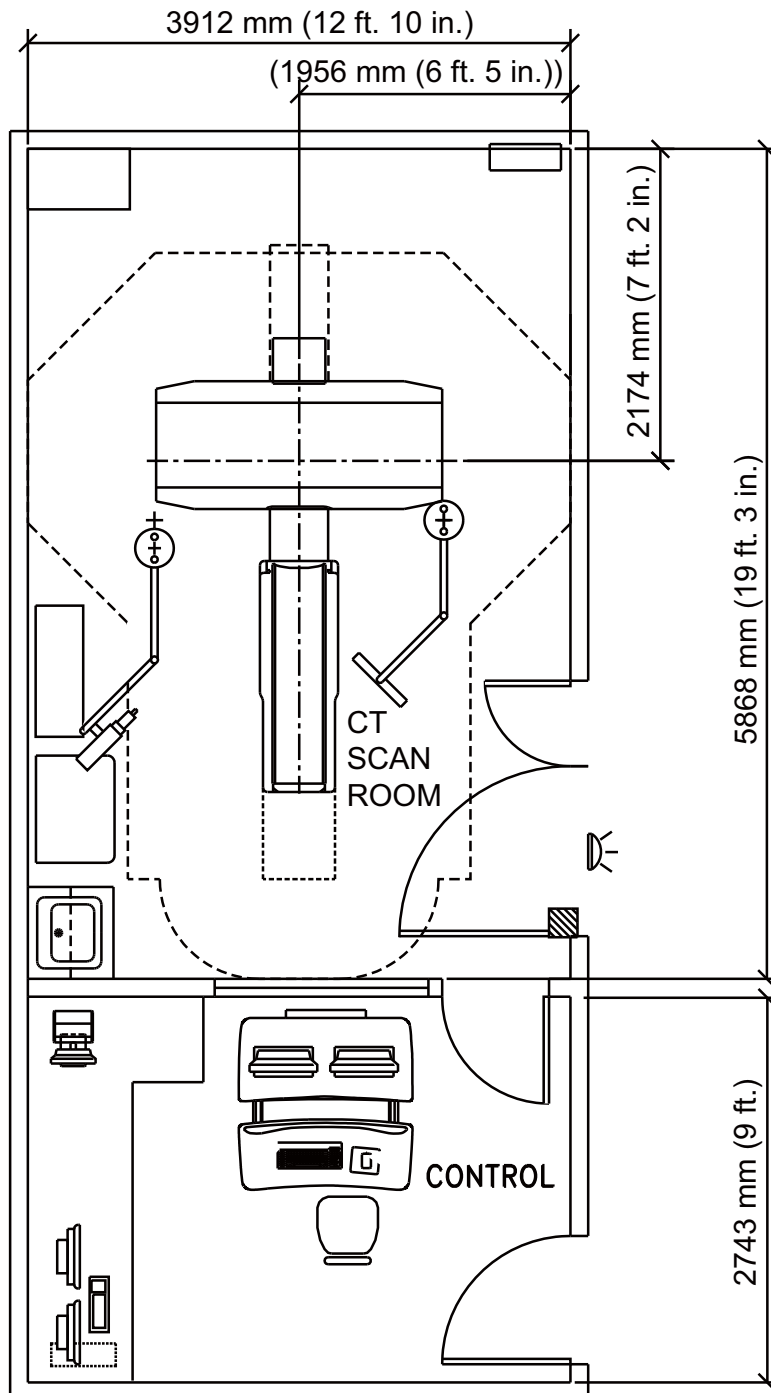


Figure E-2 Recommended Room Layout for BrightSpeed Select System


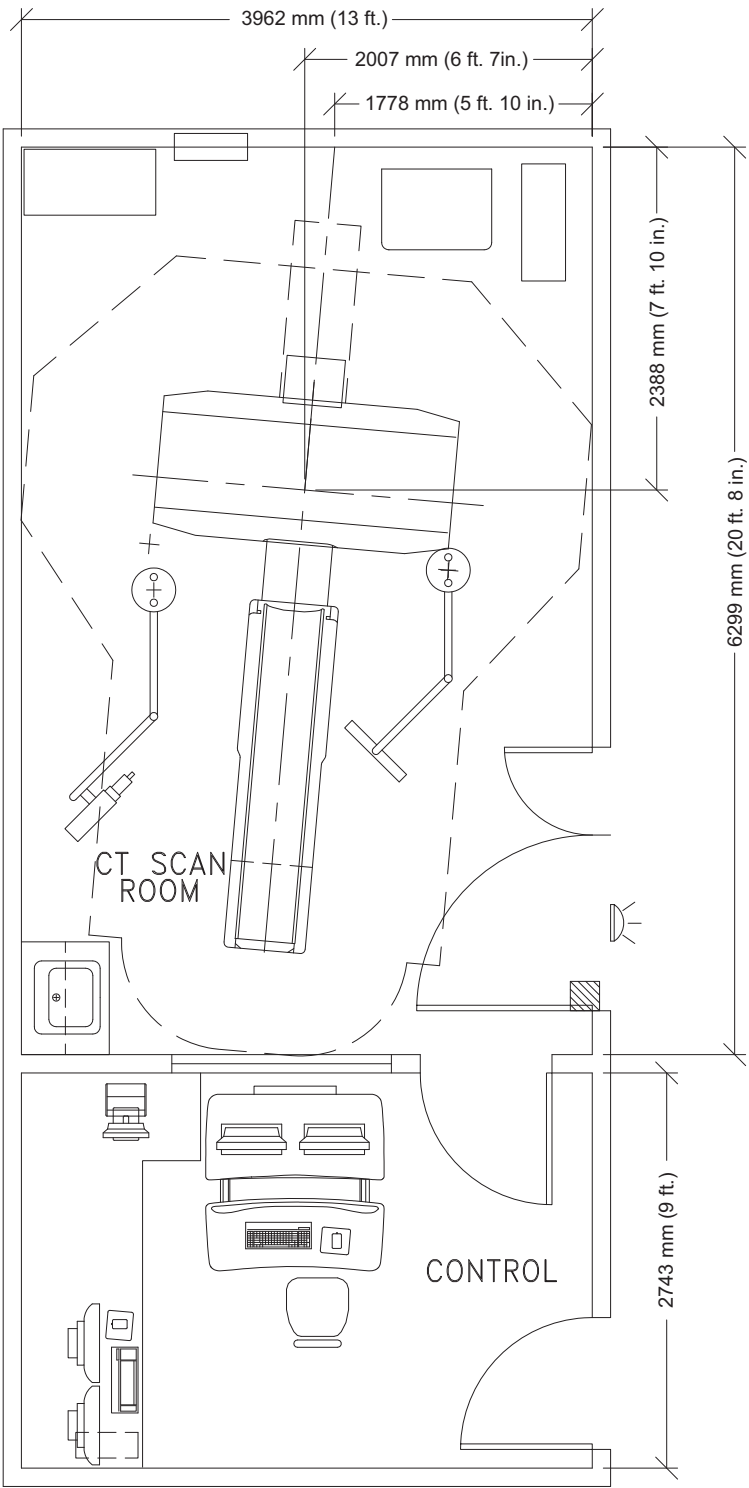


Figure E-3 Recommended Room Layout for BrightSpeed Select System



www.gehealthcare.com

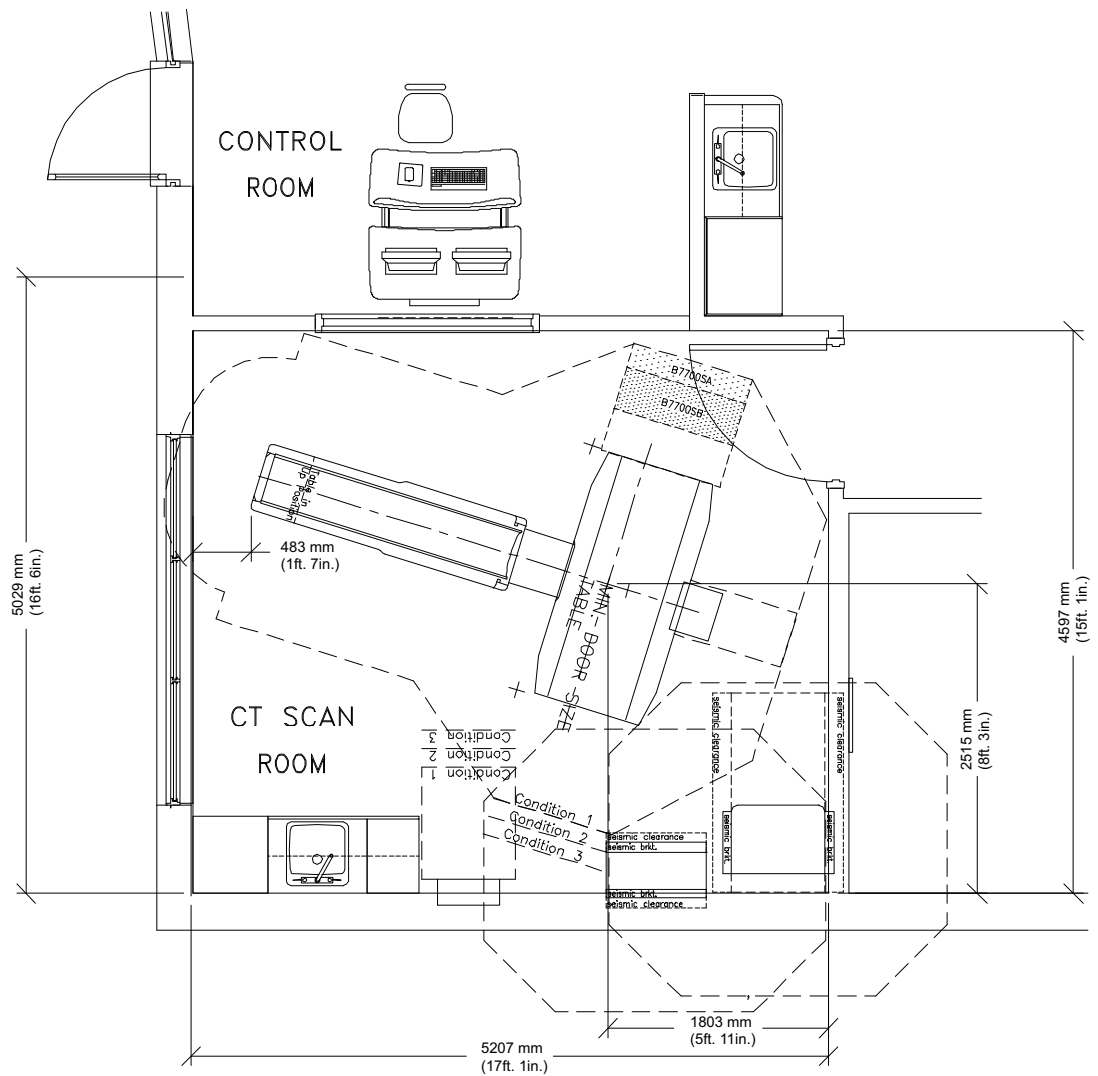