

GE Healthcare

BrightSpeed Select Series

Installation Manual

Book 2 of 2: Electrical Calibration, Integration & Testing

OPERATING DOCUMENTATION



5181248-100

Rev 7

© 2007-2008, General Electric Company

All Right Reserved

Book 2 of 2: Electrical Calibration, Integration & Testing

Pages 176 - 304

Effectivity:

The information in this manual applies to the following BrightSpeed Select Series CT Systems:

- BrightSpeed Elite Select
- BrightSpeed Edge Select
- BrightSpeed Excel Select

Table of Contents: Book 2

Chapter 4

Electrical Introduction.....	183
Section 1.0	
Introduction	183
Section 2.0	
Review Mechanical Hand Off Material.....	183
Section 3.0	
Training	183
Section 4.0	
Required FE Common Tools and Supplies.....	184
4.1 FE Calibration and Service Tool List	184
4.2 Electrical Tools	184
4.3 Image Quality Calibration Tools.....	184
4.4 Detector Service Tools	184
4.5 Optional Tools.....	185
4.6 Safety Materials	185
4.7 Cleanliness	185
Section 5.0	
Requirements/Assumptions.....	186
Section 6.0	
FE Workflow.....	186
Section 7.0	
Checklists for Completed Installation	187
7.1 System-Level	187
7.1.1 General	187
7.1.2 Optional and Regional	187
7.2 Site Clean Up.....	187
7.3 Dolly Return	187
7.4 Options	187
7.5 Paperwork (Final Activities)	188
Section 8.0	
GE and Regulatory Forms	189
8.1 All Countries	189
8.1.1 GE e-4879 Form	189
8.1.2 Product Locator Cards.....	189
8.2 U.S. Installations Only	189
8.2.1 FDA 2579 Form	189
8.2.2 System Chassis Ground Leakage Test Form	190

Chapter 5

Electrical Integration and Safety Verifications 191

Section 1.0

Electrical Power ON & Ground Checks 191

1.1	Required Tools.....	192
1.2	Initial PDU Configuration.....	192
1.2.1	Circuit Breakers.....	192
1.2.2	Relay Board	192
1.2.3	Power Switches.....	193
1.2.4	Hardware and Connection Check	193
1.2.5	Covers.....	193
1.3	Suite Emergency Off Checks	194
1.4	Line Transformer Settings.....	194
1.4.1	Requirements.....	194
1.4.2	Line Input Conditions	194
1.5	System Power Up	195
1.6	Emergency Stop Check	197

Section 2.0

Computer Integration 199

2.1	Introduction	199
2.2	System Configuration Data Sheets.....	199
2.3	Restore System State	199
2.4	Install Customer Options.....	200
2.4.1	Software	200
2.4.2	Camera	202
2.5	Shut Down Application	203
2.6	Reconfig the OC.....	203
2.6.1	Overview	203
2.6.2	Procedure.....	204
2.7	Check/Set Date and Time	210
2.8	Save System State.....	211
2.9	Applications Start-Up	211
2.10	Console Boot-up Flow Chart.....	212
2.11	Adjust Monitor	213

Section 3.0

Table Gantry Integration 215

3.1	Introduction	215
3.2	Check Alignment Lights	215
3.2.1	Room Light Adjustment.....	215
3.2.2	Turning the Alignment Lights ON	215
3.2.3	Internal Axial Lights.....	216
3.2.4	External Axial to Internal Axial Distance	216
3.2.5	Coronal Lights.....	216
3.2.6	Turn Lights OFF	216
3.3	Autovoice/Intercom Check	217
3.3.1	Requirements.....	217
3.3.2	Patient Speaker.....	217
3.3.3	Operator Console Speaker	217

3.3.4	Autovoice Volume.....	217
3.4	CT System X-Ray ON Indicators, Cautions & Warning Labels	218
3.4.1	Check And Install System Warning Labels.....	218
3.4.2	Documentation - Verification.....	218
3.5	Check Warning Labels.....	219
3.5.1	On SCIM	219
3.5.2	On Gantry	219
3.5.3	On Laser	219
3.6	Process Product Locator Cards.....	219
3.7	Verify Service Cabinet Installation	219
3.8	Check X-Ray Lights	220
3.9	Mechanical Characterization	220
3.9.1	Alignment Light Characterization.....	220
3.9.2	Table Height Characterization	220
3.10	Interference Test.....	221
3.10.1	Verify Table Elevation.....	221
3.10.2	Tests	222
3.10.2.1	Position Tilt, Move Table to Interference Limit	222
3.10.2.2	Position Table, Move Tilt to Interference Limit	223
3.10.2.3	Tilt Limits When Table Below Scan Plane Lower Limit	223

Chapter 6 Image Quality 227

Section 1.0 Introduction 227

Section 2.0 Calibration Process..... 228

2.1	Reference Procedures.....	228
2.1.1	Scanning with Service Protocols	228
2.1.2	Center Phantom.....	228
2.2	Prepare the QA Phantom	228
2.3	Calibration Process Introduction.....	228

Section 3.0 Table/Gantry Alignment Procedure 229

3.1	Time & Personnel	229
3.2	Tools and Test Equipment.....	229
3.3	Preparation	229
3.4	Procedure	229

Section 4.0 Tube Warm Up and Fast Cal..... 230

Section 5.0 Tomographic Plane Indication 232

Section 6.0 Image Series 233

6.1	Scan Protocol	233
6.2	Data Recording: Means and Standard Deviation	233

6.3	Term Definitions and Screens	233
6.4	48cm Phantom Image Series Image Performance Verification	235
6.4.1	Acquiring the 48cm Phantom Image Series	235
6.4.2	Image Performance Verification	236
6.4.2.1	Series Means (AvXo-AvXc and AvSDo) & Center Smudge	236
6.4.2.2	Band and Streak Artifact	237
6.4.3	Failure Recommended Actions	238
6.4.3.1	Series Means (AvXo-AvXc or AvSDo) Failure Recovery	238
6.4.3.2	Center Smudge Failure Recovery	239
6.4.3.3	Band or Streak Artifact Failure Recovery	239
6.5	20cm QA Phantom Image Series Image Performance Verification	239
6.5.1	Acquiring the 20cm QA Phantom Image Series	239
6.5.2	20cm QA Phantom Image Series Image Performance Verification	240
6.5.3	Failure Recommended Actions	244
6.5.3.1	1 st Image Series (4 Image MTF Average) Failure Recovery	244
6.5.3.2	2 nd Image Series (Visible Lines) Failure Recovery	244
6.5.3.3	3 rd Image Series (Visible Holes) Failure Recovery	244
6.5.3.4	4 th Image Series Failure Recovery	245

Section 7.0		
System Functional Test		247

Section 8.0		
Save System State		248

Chapter 7		
Customer Options Installation & Verification		249

Section 1.0		
CT Options		249

Section 2.0		
Install Options		250

2.1	5 inch MOD	250
2.2	Camera (Filming Device)	250
2.3	Advantage Windows Workstation (AWW)	250

Section 3.0		
DICOM Network Introduction		251

Section 4.0		
Before You Start		252

4.1	Network Physical Requirements	252
4.2	Network Identity Information	252
4.3	Scanner to DICOM Remote Hosts Network Information	252
4.4	Scanner to DICOM HIS/RIS Interface Network Information	252
4.5	Scanner to DICOM Printers Network Information	253

Section 5.0		
Declaring the System on the Hospital Network		254

5.1	Enter Configuration Routine	254
5.2	Configure Network Settings	255
5.3	Initiate System Reconfiguration	257

Section 6.0	
Declaring Remote Hosts on the CT System	258
6.1 Enter Remote Host Configuration Screen	258
6.2 Declaring Advantage NET Remote Hosts on the Scanner	258
6.3 Declaring DICOM Remote Hosts on the CT Scanner.....	259
Section 7.0	
Declaring the CT System on Remote Hosts	261
7.1 Declaring the Scanner on Advantage NET Protocol Devices/Systems	261
7.2 Declaring the Scanner on DICOM Protocol Devices/Systems	261
Section 8.0	
DICOM HIS/RIS Setup	262
8.1 Prerequisites.....	262
8.2 Loading ConnectPRO Software Option on the CT System	262
8.3 PPS Setup	265
Section 9.0	
DICOM Filming Devices Setup	266
9.1 Prerequisites.....	266
9.2 Declaring DICOM Filming Devices on the CT System	266
9.3 Troubleshooting Tips	277
Section 10.0	
Network Connections	279
Section 11.0	
Modem Setup in IIP Configuration.....	284

Chapter 8 System-Level Safety Tests 285

Section 1.0	
Patient Touch Leakage Test.....	285
1.1 Personnel Requirements	285
1.2 Overview	285
1.3 Preliminary Requirements	285
1.3.1 Tools and Test Equipment.....	285
1.3.2 Safety.....	285
1.3.3 Required Conditions	286
1.3.4 Test Conditions.....	286
1.4 Procedure	287
1.4.1 Cover Removal.....	287
1.4.2 Leakage Current Testing	287
1.5 Finalization.....	289
Section 2.0	
CT System Chassis Leakage Test	290
2.1 Personnel Requirements	290
2.2 Overview	290
2.2.1 Time.....	290
2.3 Preliminary Requirements	290
2.3.1 Tools and Test Equipment.....	290

2.3.2	Safety	291
2.3.3	Required Documentation	291
2.3.4	Required Conditions.....	291
2.4	Procedure.....	291
2.4.1	Cover Removal	291
2.4.2	Electrical Connection Removal	292
2.4.3	System Ground Wire Testing	292
2.4.4	System Power-Down.....	293
2.4.5	System Power-Up and Test	293
2.5	Finalization	293

Appendix F	
Additional Characterization Procedures.....	295

Appendix G	
System Configuration Data Sheets	297

Appendix H	
Symbols	301

Chapter 4

Electrical Introduction



NOTICE Potential for Data Loss and/or Equipment Damage

To prevent potential data loss and equipment damage, please do the following:

- When instructed, record data collected from the procedures in this chapter on GE Form e4879. For more information about this form, see [Section 8.0](#) of Chapter 4.
- Only use the Installation manual that arrives with your system for installation. Any other revisions of this manual may not exactly match your system.

Section 1.0 Introduction

Use the continuity and ground checks to verify the system power connections have not shorted to ground and that the ground and neutral connections are intact.

Section 2.0 Review Mechanical Hand Off Material

Complete the Mechanical Hand Off checklist:

- ☐ All options were installed. If not, contact your install specialist.
- ☐ Check for short ships.
- ☐ Review cable connections with mechanical team.
- ☐ Complete mechanical vendor evaluation.

Section 3.0 Training

This product requires a trained FE to proceed with the calibrations in this section.

Section 4.0

Required FE Common Tools and Supplies

4.1 FE Calibration and Service Tool List

Note: Items with a “check” (✓) are included in the Install Support Kit.
Special CT Tools used for mechanical alignments:

- | | |
|---|---|
| ✓ 5mm hex bit for 3/8" drive 6" long
(Snap-On FAML5E or equivalent) | ✓ 10mm open-end thin wrench
(Snap-On SRSM10 or equivalent) |
| ✓ 10mm hex bit for ratchet wrench (ball end)
(Snap-On FABM10E or equivalent) | ✓ 21mm open-end thin wrench
(Snap-On LTAM2124 or equivalent) |
| ✓ 2.5mm Allen hex bit for 1/4" | ✓ 14mm ball hex socket bit for 3/8" drive |
| • Standard FE Tool Kit | ✓ 14 mm hex socket bit for 1/2" drive |
| • Torque Wrench Kit | |

4.2 Electrical Tools

These tools must be calibrated yearly.

- Fluke 87 DVM or equivalent
- Clamp on amp meter

4.3 Image Quality Calibration Tools

- QA Phantom (2144715)
- IQ Cal poly phantom:
 - 35cm (2144721)
 - 48cm (2144721-2)
- Phantom Holder

4.4 Detector Service Tools

A DAS/Detector Service Kit (PN 2344539) is included in the Install Support Kit.

- Aero Duster
- AMAX Contact Cleaner

4.5 Optional Tools

Needed for warranty recalibration:

- GE HV Bleeder
- Scope with 10X probes
- Beckman CT231 clamp-on Amp probe: 46-194427P228 (Fits up to #2 size wire)
- Beckman CT232 clamp-on Amp probe: 46-194427P270 (Fits #1 size wire and larger)
- Scope Probe-to-bleeder Cable 46-219921G1 (Belden #8422 cable, 30 ft)
- Polaroid type 52 film and developer
- Radial dial indicator (mm or inches) and mounting bracket
- Caliper Dial indicators (mm or inches)
- Ground Rod Tester Clamp on Amp Probe

4.6 Safety Materials

Note: Items with “checks” (✓) are included in the Install Support Kit

- ✓ Lockout/Tagout kit, or equivalent
- ✓ Nitrile Rubber Gloves
- Safety Glasses

4.7 Cleanliness

Any dirt on the surface increases leakage current on the filter or converter cards and causes the DAS to fail the drift spec. Wear Nitrile gloves (part number 2207303-6 [large] and 2207303-7 [extra large]) when you handle the DAS, because fingerprints on the board can cause problems during humid conditions. Use only clean, new Nitrile gloves. Do not use hospital grade gloves.



NOTICE NEVER USE AN ERASER TO CLEAN ANY PART OF THE DAS.
DUST COVERS REQUIRED for installations on construction sites.

Section 5.0

Requirements/Assumptions

- This procedure will be performed by an appropriately trained engineer.
- All stations in a suite must have the same SUITE ID.
- You need the name of all hosts in the suite and their corresponding Internet/Ethernet numbers.
- You need the Internet (IP) addresses the first time you execute a reconfig on the system.
 - When you connect the system to a network, contact the system administrator to obtain the IP addresses for all the computers in the suite.
 - You also need an IP address for each gateway (second) ethernet board in any OC or IC.
 - You can use the default internet number on stand-alone systems (stand-alone = not connected to any network).
- For more detailed information and instructions regarding Network Integration, see [Chapter 7 on page 249](#).

Section 6.0

FE Workflow

- 1.) Review mechanical hand-off material.
- 2.) Obtain required FE common tools and supplies.
- 3.) Perform electrical power-on and ground checks.
- 4.) Gather all customer information needed for reconfiguration.
- 5.) Perform computer integration.
- 6.) Complete Table/Gantry integration.
- 7.) Complete the calibration process.
- 8.) Perform the Table/Gantry Alignment procedure.
- 9.) Perform tube warm-up and fast calibration.
- 10.) Complete tomographic plane indication.
- 11.) Run image series tests.
- 12.) Run system functional test.
- 13.) Create system state DVD.
- 14.) Perform the Patient Touch Leakage Test.
- 15.) Perform the CT System Chassis Leakage Test, as required by local code.
- 16.) Complete installation and verification of any customer options.
- 17.) Complete and return GE Form e4879 Installation Data Verification, for all installations.

Section 7.0

Checklists for Completed Installation

Complete the installation tasks listed below and check the appropriate boxes here and on the GE e-4879 form to verify the completion of these tasks. Section 8.1 contains an explanation of the GE e-4879 form, which the FE must complete and submit for ALL installations.

7.1 System-Level

7.1.1 General

- ☐ HVAC system is operational and environmental data reported on the GE e-4879 form.
- ☐ System realignments completed, if required.
- ☐ Broadband installed and operational.
- ☐ Power and ground audit completed.

7.1.2 Optional and Regional

- ☐ Seismic mounting kit installed, if required in your area.
- ☐ Generator recalibration completed, if necessary.
- ☐ Collimator recalibration completed, if necessary.

7.2 Site Clean Up

- ☐ All DVDs for customer options placed in the GE service cabinet.
- ☐ All system software and service tools placed in the GE service cabinet.
- ☐ System cleaned and nicks touched-up with paint.
- ☐ Installation site cleaned and all trash properly disposed.

7.3 Dolly Return

- ☐ Return of dollies arranged and dolly pick-up confirmed.

7.4 Options

Check the appropriate boxes here and on the GE e-4879 form to verify the installation and proper functionality of all customer-ordered options.

- ☐ Injector installed and operational.
- ☐ Advantage Windows Workstation installed and functional tests completed.
- ☐ Advantage 4D installed and functional tests completed.
- ☐ Filming/Camera/DASM installed and operational.
- ☐ Modem installed and functional tests completed.
- ☐ UPS installed and functional tests completed.

- ☐ Network items installed and functional tests completed.
- ☐ Customer software options installed and operational.
- ☐ Teleradiology connections completed. See Section 10.0 of Chapter 7.
- ☐ Remote monitor installed and operational.
- ☐ Bar code reader installed and operational.
- ☐ Cardiac monitor and stand installed and operational.

7.5 Paperwork (Final Activities)

- ☐ GE e-4879 form completed; see Section 8.1. (Required for installations in ALL countries.)
- ☐ FDA 2579 form completed; see Section 8.2. (Required ONLY for U.S. installations.)
- ☐ Any PQRs or PSRs that you encountered have been reported.
- ☐ All FMIs for system completed, if necessary.
- ☐ All dispatching activities (03-04-10 codes) completed.
- ☐ Customer acceptance checks completed.
- ☐ System transfer completed and appropriate GE Healthcare personnel notified.
- ☐ All outstanding customer installation issues have been addressed.

Section 8.0

GE and Regulatory Forms

Field Engineers must complete and submit the documents listed in Section 8.1 for ALL installations, regardless of the country. In addition, for installations performed within the United States, Field Engineers must ALSO complete and submit the documents listed in Section 8.2.

8.1 All Countries

8.1.1 GE e-4879 Form

The Field Engineer should:

- 1.) Locate the GE e-4879 form on the Service CD.
- 2.) Complete the form.
- 3.) E-mail the completed form to the following:
 - a.) HHS Administrator

8.1.2 Product Locator Cards

The Field Engineer should:

- 1.) Enter the Product Locator Card information on the Product Locator Web site. Go to the following address to access the site: http://gib.gehealthcare.com/gib/gib_entry.jsp
- 2.) Leave ONE (1) Product Locator Card (or a copy) at the customer site for EACH piece of equipment installed there.

Note: CT Manufacturing completes GE HHS Data Sheets and provides them to the HHS Administrator.

8.2 U.S. Installations Only

8.2.1 FDA 2579 Form

The Field Engineer should:

- 1.) Download the FDA 2579 form from the HHS Support Central Web site:
http://supportcentral.ge.com/products/sup_products.asp?prod_id=16442
- 2.) Complete the form.
- 3.) E-mail the completed form to the HHS Administrator.

Note: Do NOT print this form after completion. The HHS Administrator will e-mail a printable version to the FE for customer site records.



NOTICE

Some states require a State Registration Number to complete this form. For any questions concerning your state, contact the HHS Administrator or check the HHS Support Central Web site. Some states may also require additional information and test information. For instructions, contact the Project Manager of Installation.

8.2.2 System Chassis Ground Leakage Test Form

The Field Engineer should:

- 1.) Locate the System Chassis Ground Leakage Test form on the Service CD.
- 2.) Complete the form, if required in your area.
- 3.) Forward the results as instructed on the form.

Chapter 5

Electrical Integration and Safety Verifications



NOTICE Potential for Data Loss and/or Equipment Damage

To prevent potential data loss, please do the following:

- When instructed, record data collected from the procedures in this chapter on GE Form e4879. For more information about this form, see [Section 8.0](#) of Chapter 4.
- Only use the Installation manual that arrives with your system for installation. Any other revisions of this manual may not exactly match your system.

Section 1.0

Electrical Power ON & Ground Checks



WARNING



THIS PROCEDURE MEASURES POTENTIALLY HAZARDOUS VOLTAGES. USE AND FOLLOW LOCKOUT/TAGOUT PROCEDURES.

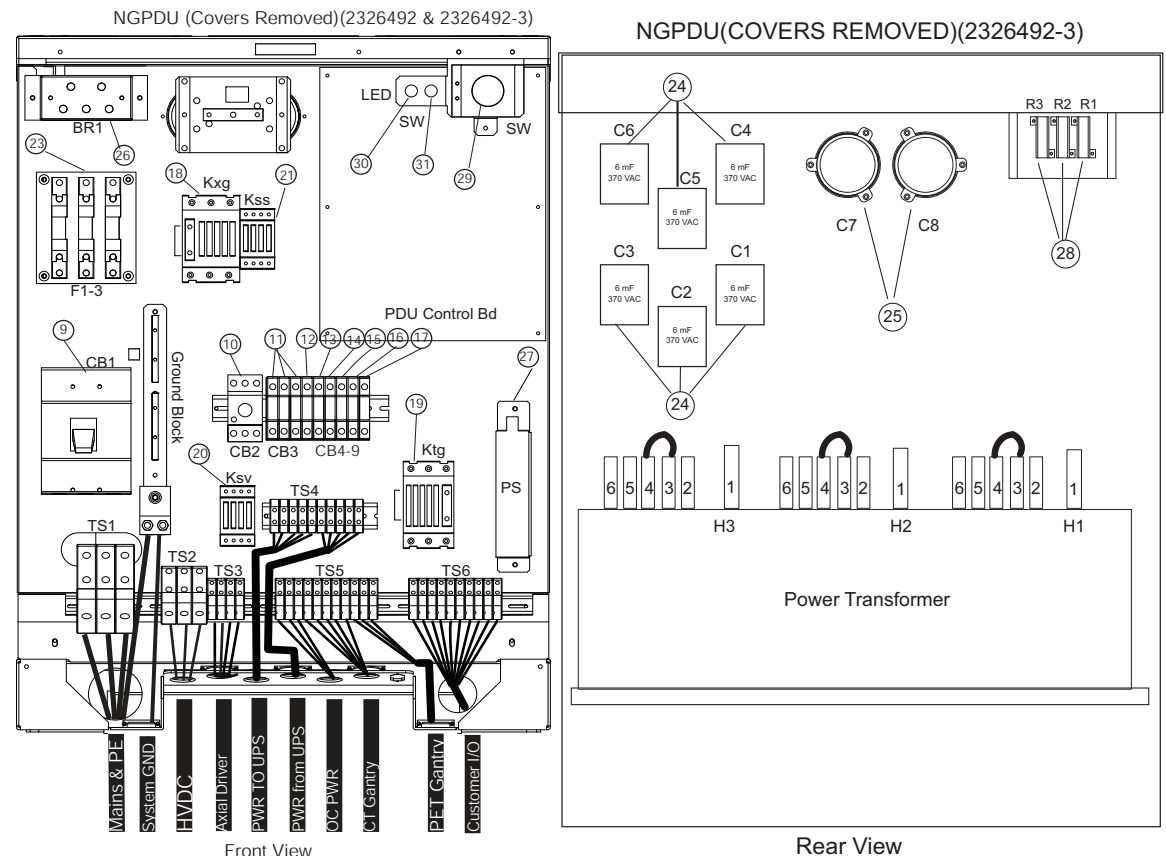


Figure 5-1 NGPDU

1.1 Required Tools

- Multimeter with a rating of at least 1000 volts
- Multimeter leads with a rating of at least 1000 volts

1.2 Initial PDU Configuration



WARNING



THIS PROCEDURE MEASURES POTENTIALLY HAZARDOUS VOLTAGES. USE AND FOLLOW LOCKOUT/TAGOUT PROCEDURES.

1.2.1 Circuit Breakers

Set all circuit breakers to OFF

1.2.2 Relay Board

- 1.) Set SW to the normal position.
- 2.) When system is already, three lamps are both lighting (refer to [Figure 5-2](#)).

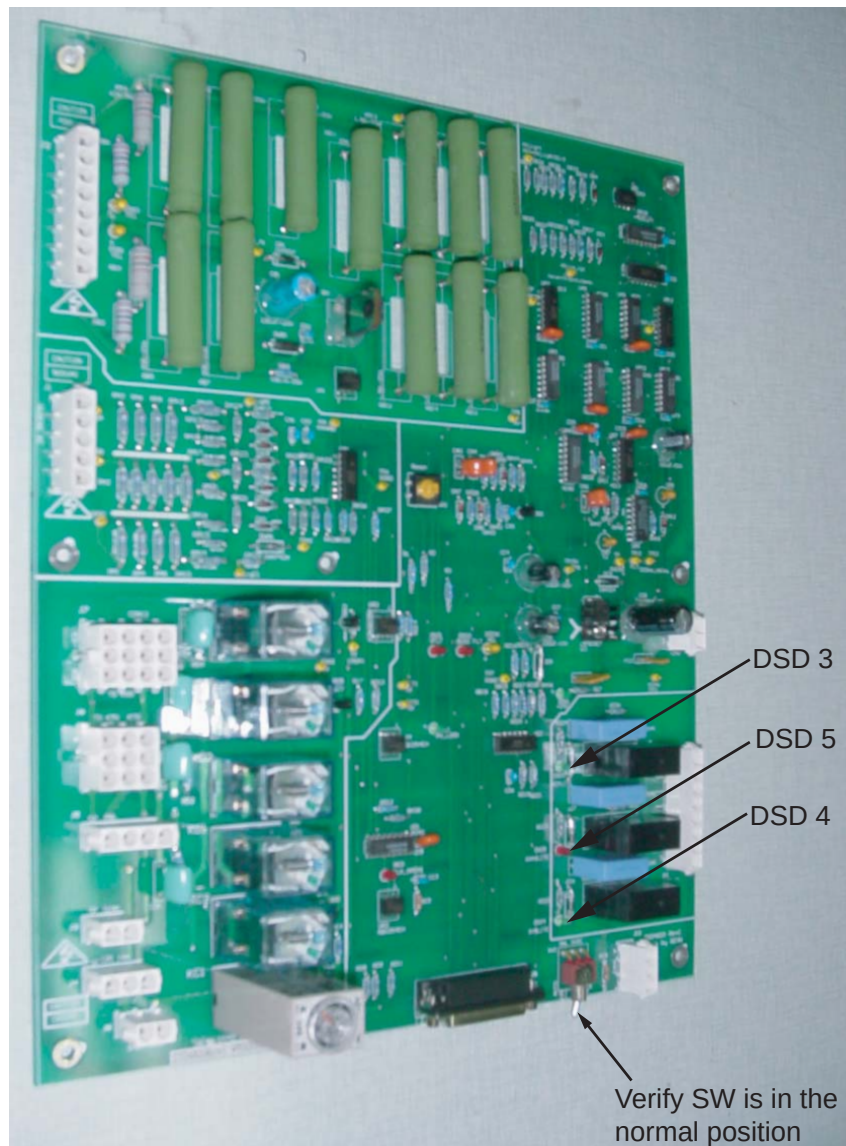


Figure 5-2 NGPDU Control Board

1.2.3 Power Switches

Turn OFF all power switches on all subsystems

- Gantry
- DAS
- Power Pan Breaker
- Table
- Console

1.2.4 Hardware and Connection Check

Use this step to check mechanical connections and tighten anything that may have shaken loose during shipment. Verify **all** hardware and connections in the PDU are securely fastened.

1.2.5 Covers

Install, or verify the presence of, all the lexan safety covers.

1.3 Suite Emergency Off Checks



WARNING

VERIFY ALL PERSONNEL HAVE CLEARED THE SYSTEM BEFORE YOU TURN ON WALL POWER.

- 1.) Turn wall power ON to the PDU.
- 2.) **Press the suite emergency off button and verify it turns off wall power to the PDU.**
[Typically, this red palm button is located on the wall close to the console, within the scan suite.]
- 3.) Verify that all "Emergency Off" button are working properly.
- 4.) Leave power "OFF"

1.4 Line Transformer Settings

1.4.1 Requirements

- 1.) The PDU is shipped configured for 480VAC.
- 2.) Complete only if your site uses a voltage other than 480VAC.
- 3.) If PDU is configured for 480VAC, go to [1.5](#). Otherwise, proceed to [Section 1.4.2](#).



WARNING

MAKE SURE YOU TURNED OFF, TAGGED AND LOCKED THE MAIN WALL POWER BEFORE YOU CHANGE TAPS. FAILURE TO DISCONNECT POWER AT MAIN INPUT MAY RESULT IN ELECTROCUTION. TURN OFF WALL POWER TO CONNECT OR MOVE METER LEADS, OR TO REMOVE OR INSTALL COVERS.

Taps should be shipped as shown for 480 VAC only. For all others, you must move the taps. The tap check should be completed by the mechanical installer.

1.4.2 Line Input Conditions

- 1.) Monitor the No Load Line to Line Voltage at L1, L2, L3, during the workday. Do not record this data during "brown out" conditions.
- 2.) After you determine the nearest nominal line, verify the tap connections match (refer to [Table 5-1](#) and [Figure 5-3](#) for tap locations).

Transformer Taps and Jumpers

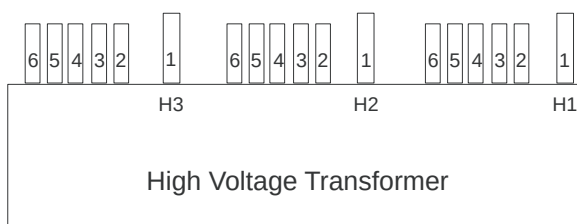


Figure 5-3 PDU Tap Positions (Rear)

- 3.) Verify that the No Load Line to Line Voltage never falls outside the corresponding minimum and maximum values listed in [Table 5-1](#).
- 4.) Use a 0-750 AC voltmeter of 3/4% accuracy to measure the line-to-line voltages at L1, L2, and L3.
 - Verify the highest line-to-line voltage does not exceed 1.02 times the lowest voltage.
 - **Example:** If the lowest voltage equals 474, the highest voltage should not exceed $474 \times 1.02 = 483.5$ volts.



THIS PROCEDURE MEASURES POTENTIALLY HAZARDOUS VOLTAGES. USE AND FOLLOW LOCKOUT/TAGOUT PROCEDURES.

NO LOAD Line to Line Voltages		TAP CONNECTIONS (All 3 phases must have same the configuration)		
Nominal	Maximum Range (10%)	Phase A Connection	Phase B Connection	Phase C Connection
480V*	432 to 528*	3-4*	3-4*	3-4*
460V	414 to 506	3-5	3-5	3-5
440V	396 to 484	3-6	3-6	3-6
420V	378 to 462	2-4	2-4	2-4
400V	360 to 440	2-5	2-5	2-5
380V	342 to 418	2-6	2-6	2-6
240V**	216 to 264**	1-4**	1-4**	1-4**
220V**	198 to 242**	1-5**	1-5**	1-5**
200V**	180 to 220**	1-6**	1-6**	1-6**
* Factory Default				
** 2326492-3 PDU only				

Table 5-1 PDU Line Tap Connections

1.5 System Power Up



CAUTION Verify all personnel have cleared the system before you turn on wall power.

- 1.) Turn ON A1 breaker panel.
- 2.) Turn ON all system power switches and breakers (table, gantry, DAS, PDU, console).

SUB-SYSTEM POWER-UP

- 1.) Turn ON switch S3 in the table (120vac 24hr power).
- 2.) Turn ON the Gantry **120VAC**. (Light should turn on.)
- 3.) Turn **AXIAL DRIVE ENABLE** ON. (Light should turn on.)
- 4.) Turn **HV DC ENABLE** ON. (Light should turn on.)
- 5.) Push the Service Switch Panel reset button (see [Figure 5-4](#)).

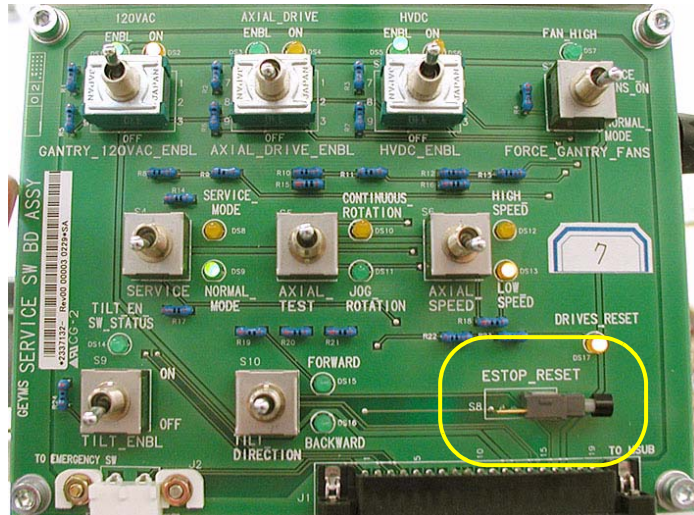


Figure 5-4 Service Switch Panel

AXIAL ENABLE SWITCH TEST

- Note:
- 1.) Turn OFF axial drive enable switch AXIAL DRIVE on the Service Switch Panel.
 - 2.) Clear the gantry area for rotation.
 - 3.) Press the alignment light push button.
 - 4.) Verify that the gantry did not rotate.
 - 5.) Turn ON axial drive enable switch AXIAL DRIVE on the Service Switch Panel.

ROTATION SAFETY CHECKLIST

- 1.) Turn OFF axial drive enable switch AXIAL DRIVE.
- 2.) Turn OFF HVDC enable switch.
- 3.) Press red E-STOP button.
- 4.) Manually rotate the gantry 360 degrees. (Keep one finger on E-stop button.)
 - Listen for any interference between the rotating and stationary parts. (Correct any interference problems.)
 - Listen for any loose parts. (Tighten parts as needed.)
- 5.) Turn ON axial drive enable switch AXIAL DRIVE.

WARNING MAKE SURE THERE ARE NO OBSTRUCTIONS AROUND THE GANTRY. PRESSING THE ALIGNMENT LIGHT PUSHBUTTON WILL CAUSE THE GANTRY TO ROTATE.

- 6.) Press the alignment light push button.
- 7.) Verify that the gantry rotates.
- 8.) Perform a 4 second X-ray OFF scan.

NOTICE During the scan, it may be necessary to enter the scan room, to obtain a better listening position. If so, keep a finger on one of the four E-STOP buttons (on the gantry), to quickly stop the gantry, if necessary.



- a.) From the console, click on the SERVICE DESKTOP icon.
- b.) Select DIAGNOSTICS.
- c.) Select DIAGNOSTIC DATA COLLECTION
- d.) Set the scan time to 4.00 seconds
- e.) Leave the door open. (This makes it easier to hear any loose or interfering parts.)
 - * Listen for any interference between the rotating and stationary parts.
(Correct any interference problems.)
 - * Listen for any loose parts.
(Tighten parts as needed.)
- 9.) After completing the 4 second scan, repeat steps 9a through 9e, with the following scan times:
 - 2.0 second scans
 - 1.0 second scans
 - 0.7 second scans
 - 0.5 second scans

1.6 Emergency Stop Check

- 1.) Use the gantry push-buttons to advance the cradle about 0.5m (2ft) from the home position.
- 2.) Press one of the E-STOP buttons on the gantry.
- 3.) Make sure the TABLE POWER shuts off, and the green LED flashes.
- 4.) Depress one of the table elevation buttons, to verify the emergency stop disabled table elevation.
- 5.) Depress one of the cradle drive buttons, to verify the emergency stop disabled the cradle drive.
- 6.) Press one of the RESET buttons to turn on X-RAY DRIVES POWER. (120 VAC LED stops flashing.)
- 7.) Press the other E-STOP button on the gantry.
 - a.) Make sure the Table Power shuts off.
 - b.) Manually move the cradle to the home position to make sure the cradle released.
 - c.) Make sure the cradle latches securely in the home position.
- 8.) Press one of the RESET buttons to turn on X-RAY DRIVES POWER.
- 9.) Press one of the four table tape switches to make sure the table down motion stops. Repeat with the three remaining table tape switches.
- 10.) Press the console emergency stop switch; make sure the Table Power shuts off.
- 11.) Press one of the RESET buttons to turn on X-RAY DRIVES POWER. (See [Figure 5-4](#)).



Figure 5-5 Reset buttons on Gantry and Service Switch bank.

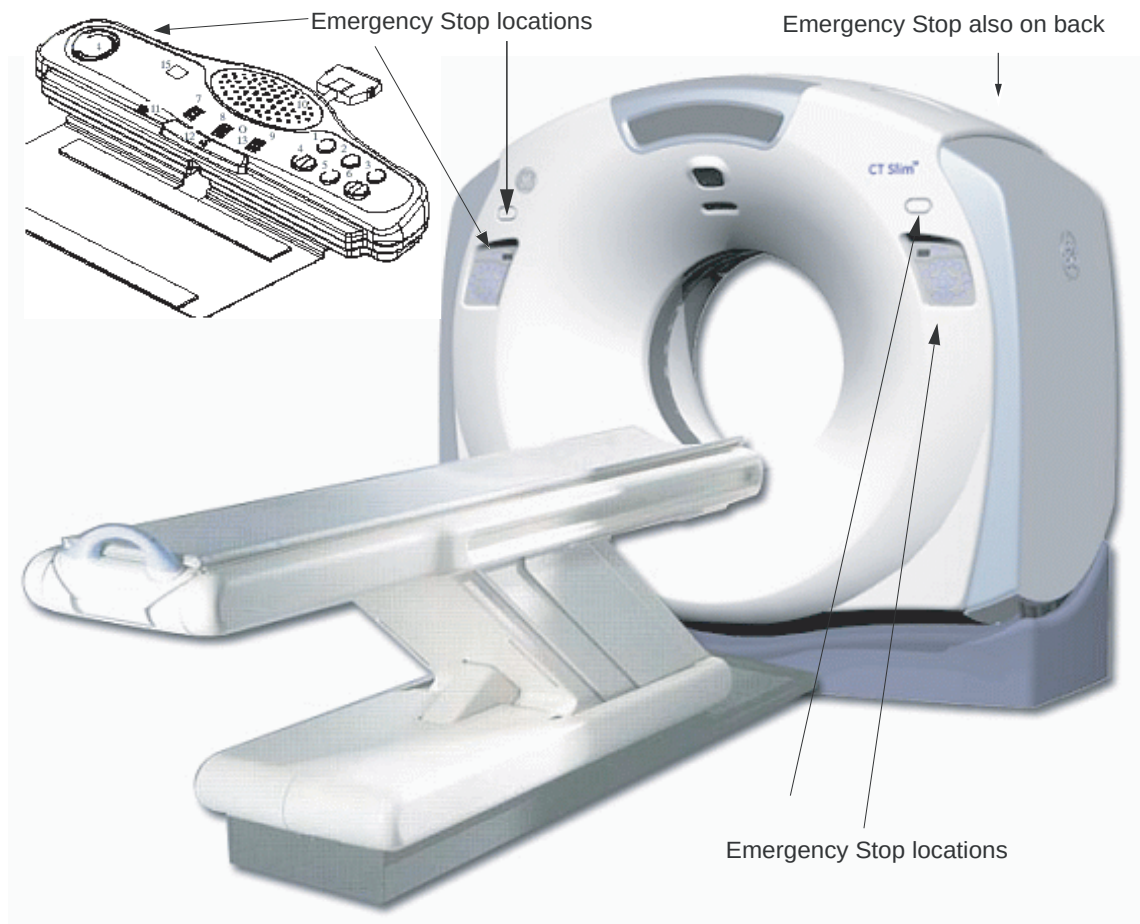


Figure 5-6 Emergency Stop button positions

Note: Emergency Stop buttons are located on the front and rear of the gantry (8 in all). They are also located on both sides of the table base (4 in all) as noted in [Figure 5-6](#).

Section 2.0 Computer Integration

2.1 Introduction

This Section describes the reconfiguration, system state restore, options, and monitor adjustment procedures.

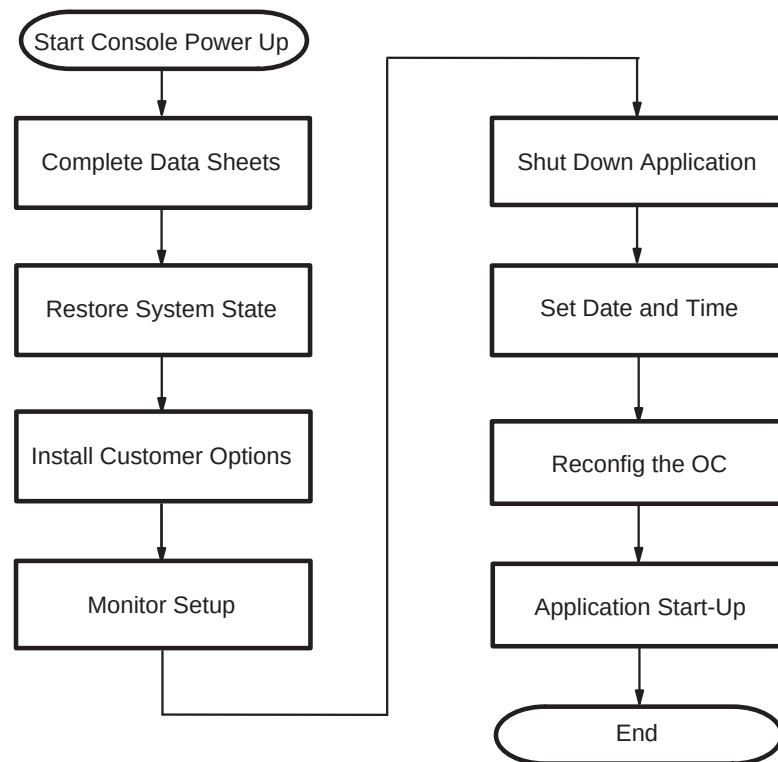


Figure 5-7 Computer Integration Process Overview

2.2 System Configuration Data Sheets

For convenient removal and use during installation, System Configuration Data Sheets appear in [Appendix G](#). Please locate and complete to them at this point during installation.

2.3 Restore System State

Tools Required: None

Your system should have a system state DVD, located in the Gantry accessories' box. The system state DVD contains:

- Collimator Characterization
- Phantom Calibrations
- Gen Cal
- Other Data

The installation process uses all the system state files. At this time, use the system state DVD to restore all files.

If you cannot locate an existing system state DVD, you must recalibrate your system.



- 1.) If you are not on the Service Desktop, click on the SERVICE DESKTOP icon.



- 2.) Click on the UTILITIES icon.
- 3.) Select SYSTEM STATE.
- 4.) Insert the DVD into the DVD drive.
- 5.) Select CHARACT.
- 6.) Select CALS.
- 7.) Select RESTORE to restore the system characterization and phantom calibration files to the system.

Note: **Restore State can take as long as ten minutes, although typical times average about three minutes. When Restore State completes, dismiss the tool, and proceed to the next section.**

If any error should occur during the restore process, see the Software Load Procedure manual (Load From Cold) for information regarding possible error messages and their recovery.

- 8.) Click NO for Reset Scan Hardware popup message.
- 9.) Select DISMISS.

2.4 Install Customer Options

Tools Required: None

2.4.1 Software

Note: **Your system may have a DVD that contains customer purchased options. If your system has an options DVD, install it at this time--otherwise skip this section.**

Ensure that the options DVD is NOT write protected at this time. The initial install requires that the DVD be write enabled; subsequent installs can be done with the DVD write protected.



- 1.) If you are not on the Service Desktop, click on the SERVICE DESKTOP icon.



- 2.) Click on the CONFIGURATION icon.
- 3.) Select INSTALL OPTIONS.

4.) Select **INSTALL OPTIONS**.

An Options Window appears (Figure 5-8):

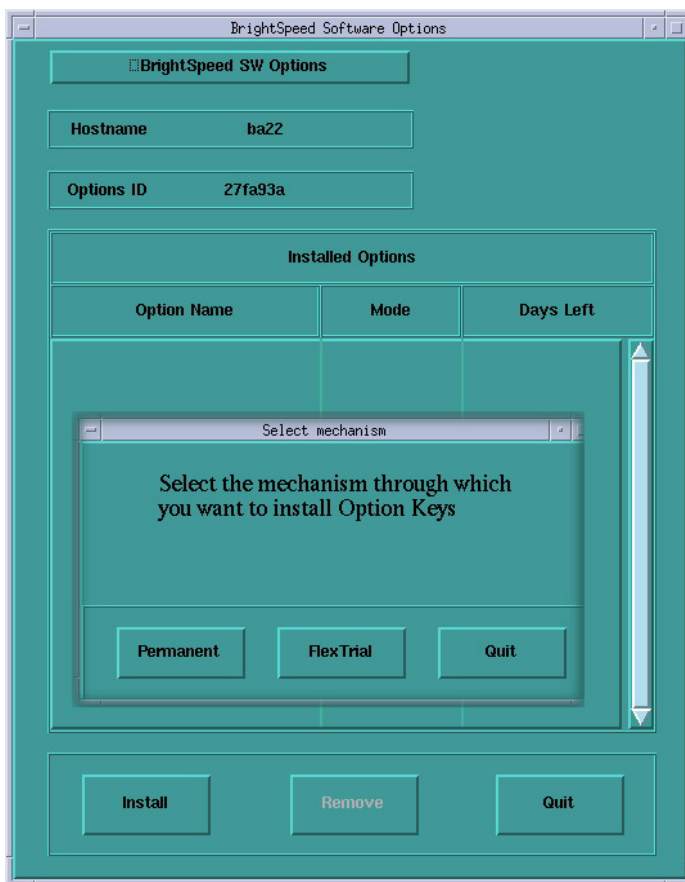


Figure 5-8 Options Window when First Selected

- Check the FDO to see what options were ordered.
- Compare FDO options to those on the Option DVD.
- If different contact your local sales representative.

- 5.) Insert the options DVD into the DVD drive and click on OK. (If you do not have an options DVD, click on OK anyway, wait for the abort pop up, then abort the process.)

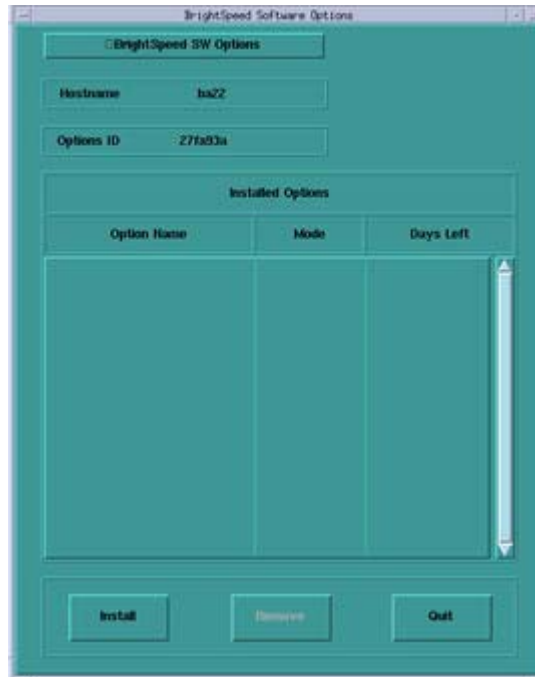


Figure 5-9 Example: Options Window

- 6.) Select all of the options in the left-hand column to install the corresponding software.
- 7.) Select INSTALL. A box may appear while the options are loading. When an option is displayed in the Installed Options list, then installation of that option is complete. Note that some options take a fraction of a second to install, while options like 3D may take a half minute (due to the fact that they are installing software).
- 8.) After the options are installed, select QUIT.
- 9.) Select OK.
- 10.) Remove the DVD and write protect the side with options.
- 11.) When the system prompts to Reboot, click YES, and reboot the system to complete the installation.

2.4.2 Camera

Tools Required:

- Small flat blade screw driver
- Data collected from data sheets ([See Camera Application Configuration on page 298.](#))
- Software Load Procedures manual
- System Service manual.

Note: If a DASM is required, the DASM hardware must be installed before proceeding.
For details on camera configuration, refer to the Software Load Procedures manual.
For details on troubleshooting the camera, refer to the System Service manual.



- 1.) Click on the SERVICE DESKTOP icon.



- 2.) Select CONFIGURATION icon.

- 3.) Select INSTALL CAMERA.
- 4.) Select ADD.
- 5.) Select DASM or DICOM.
- 6.) Follow procedures on the screen.
- 7.) Return to Home Page.

2.5 Shut Down Application

Tools Required: None

Standard Level

If Applications is currently running, you must shutdown system applications.

- 1.) Click on the SERVICE DESKTOP icon. 
- 2.) On the desktop toolbar select UTILITIES icon. 
- 3.) Select APPLICATIONS SHUTDOWN (to bring down applications only).

Super User Level

- 1.) Open a UNIX Shell window.
- 2.) Type: **su** - ENTER at the prompt.
- 3.) Type the root (super user) password: **#bigguy**

2.6 Reconfig the OC

Tools Required: None

Note: The document collector box that arrived with your system contains the *Software Installation Procedures* manual, which documents the reconfiguration procedure in more detail.

2.6.1 Overview

On the following screens, you should make the changes necessary, pressing the corresponding button at the top of the screen to move from screen to screen. When you are done, you can either press the ACCEPT button to start the reconfiguration process, or press the QUIT button to exit without changing the system configuration.

While the reconfiguration is going on, messages are displayed in a shell window that closes when reconfiguration is complete. Should you later want to review the reconfiguration output, it is logged to the following file:

`/var/adm/install.log.YYYYMMDDWWHHMMSS`

Where

YYYYMMDDWWHHMMSS is the Date/Time that the reconfiguration was started.

To view the file, type: **`more /var/adm/install.log.YYYYMMDDWWHHMMSS`**

It is possible to abort the reconfiguration while entering information on the reconfiguration screens. Simply press the QUIT button at the top of the screen. There is NO safe way to abort the reconfiguration after pressing the ACCEPT button. If the entries made in the screens were incorrect, DO NOT try to stop the reconfiguration, instead wait for it to complete, and rerun reconfig, entering the correct parameters.

2.6.2 Procedure

- 1.) Change directory to scripts:
Type: **cd /usr/g/scripts** ENTER at the prompt.
- 2.) Launch the Install utility:
Type: **reconfig** ENTER at the prompt.
The OC displays the Install Utility Window as shown in [Figure 5-10](#).

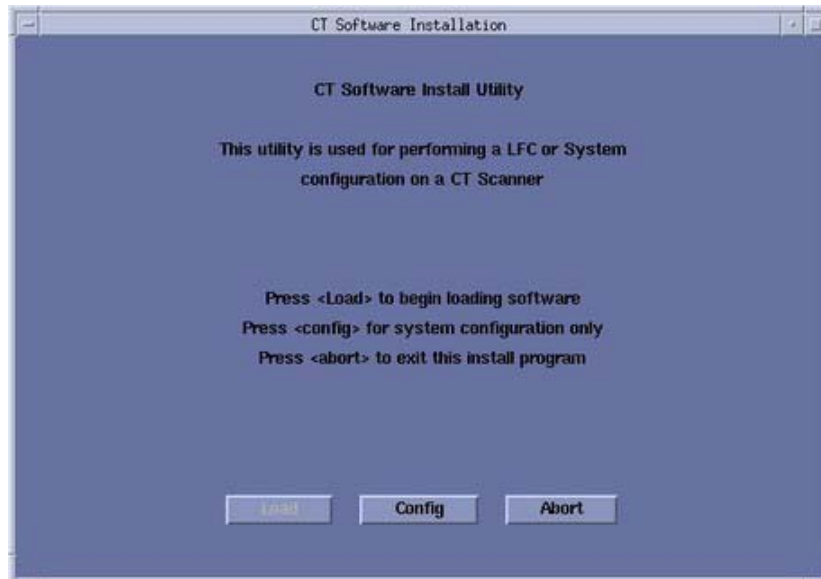


Figure 5-10 Install Utility Window

- 3.) Click on the CONFIG button.
The OC displays the System Configuration - System Settings Screen as shown in [Figure 5-11](#).

Comment: The following pages show the screens that are used to change the configuration of the system. These screens are the same as those used for the Software Configuration during Load From Cold. The actual screens will vary depending on the current configuration of your system.

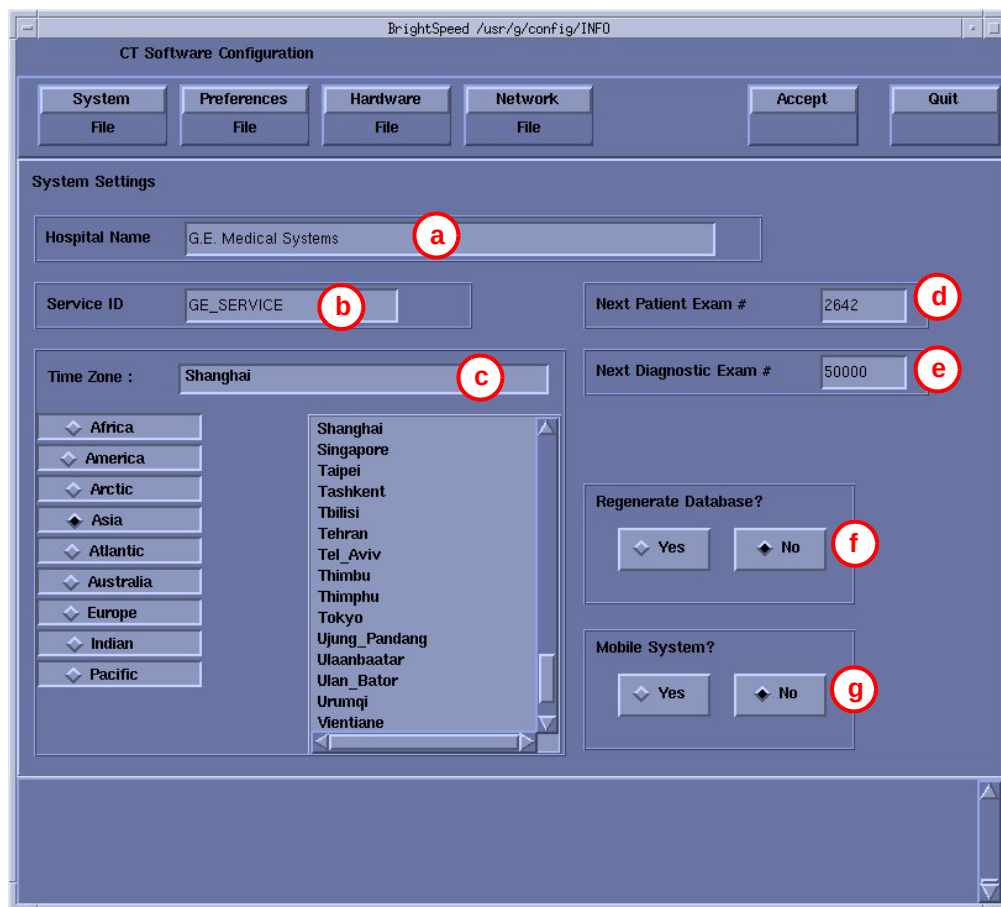


Figure 5-11 System Settings Screen

4.) Configure System Settings

- a.) Hospital Name (Figure 5-11, item 1) configures the name that will show up on images produced by this scanner. *Example:* ST MARYS HOSPITAL
- b.) Service ID (Figure 5-11, item 2) is issued by the Service organization.
Example: 262785CT2 (no spaces)
- c.) Select the Time Zone for the site.

Note: Use the scrollbar at the bottom of the time-zone selection list to view the entire description of the time-zone you are about to select, to ensure that you are selecting the correct time-zone for your location.

If the time-zone of your location is not in the list above, select one of the universal times in the selection menu. In this case, automatic changes for daylight savings time will not take effect. See Load from Cold manual, if you require more information regarding time-zone setting & selection.

- d.) Next Patient Exam # configures the next Exam number the scan user interface will use. At initial system installation, type: **1**
 - e.) Next Diagnostic Exam # *Customer Selected.*
 - f.) Regenerate Database
 - g.) Mobile System Select to tell the software if this CT is in a mobile environment or not.
- 5.) Select the PREFERENCES button to display the Preference Settings Screen as shown in Figure 5-12.
- 6.) Regenerate Database (refer to Figure 5-12) determines whether the Scan Database will be regenerated during reconfiguration. You should select NO.

Note: If YES is selected, then all scan data will be deleted on the Scan Data Disk. To avoid accidental regen, NO is selected by default each time reconfig is run. Regenerating the Scan Database destroys all the scan and calibration files on the Scan Data Disk. Make sure that all images have been reconstructed before regenerating the Scan Database.

This procedure is normally only needed when a Scan Data Disk is replaced or reformatted. See Load from Cold manual, if you require more information regarding Scan Database Regeneration.

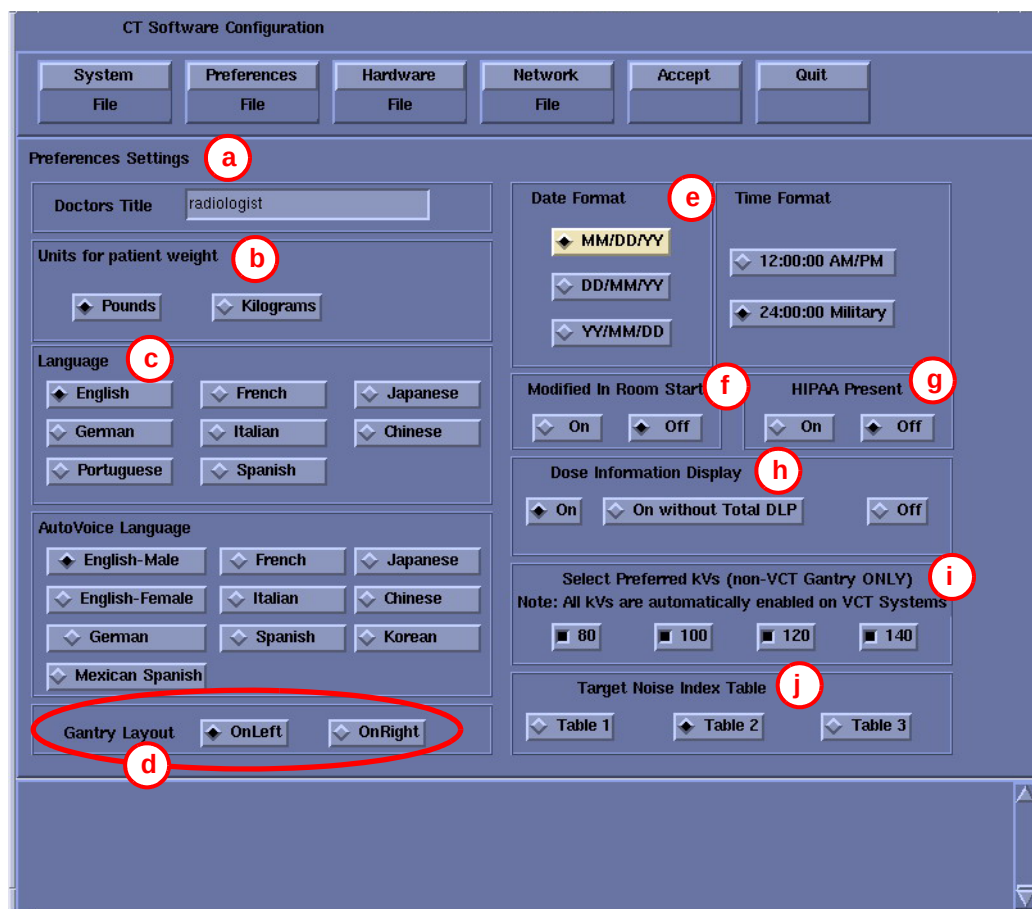


Figure 5-12 Preferences Setup Screen

7.) Configure Preferences Settings

- a.) Doctors Title Enter the title for the Doctor. (eg. radiologist)
- b.) Units for Patient Weight Tells the software whether pounds or kilograms are being used.

Note: Use Direction 5221102-1EN (found in the keyboard collector kit shipped with the system) to complete the language selection.

- c.) Language configures the language to be displayed on the Application screens.
 - 1.) Review the language matrix (in direction 5221102-1EN) and identify the appropriate language for your country.
 - 2.) Ask the radiology manager, or equivalent, to agree to the user interface (UI) language and keyboard requirements listed for that country (per 5221102-1EN), OR choose another language for the UI and keyboard.
 - 3.) Record the information on GE form e4879.

Note: d.) Select Gantry Layout according to gantry/table layout in the hospital.
Select ONLEFT if gantry is laid to the left of table (viewed from the control room), otherwise, select ONRIGHT.

- e.) Date Format configure the format in which the date will be displayed on the images. Time Format configure the format in which the time will be displayed on the images
- f.) Modified In Room Start: Select ON for Japan sites, OFF for other sites.

- g.) Set HIPAA Present to [OFF] unless the customer requests HIPAA to be on.
 - h.) Select the site preferred Dose Information Display option for the site to use in monitoring calculated Patient Dose:
 - * Select ON (full CTDiw Display)
 - * Select ON WITHOUT TOTAL DLP (no Dose Length Product Display), or
 - * Select OFF (no CTDiw Display)
 - i.) Preferred Fast Cal KV configures the preferred kV that the Fast Cal Routine will calibrate (80, 100, 120, 140 in the Selected Preferred Fast Cal KV field). The default selections are 120 and 140.
 - j.) Select Target Noise Index Table: default is [TABLE 2].
- 8.) Select the HARDWARE button to display the Hardware Settings Screen (Figure 5-13).

Comment:

These kVs should include all kVs that the site uses for patient scanning. Deselecting All Preferred FastCal KVs is the same as selecting ALL the Preferred FastCal KVs



Figure 5-13 Hardware Settings Screen (example only - actual screen may vary)

- 9.) Configure Hardware Settings
Verify correct hardware type is shown. If not correct, select correct DAS Type (**4/8 Slice** only), Table Type and Numbers of IGs.
16 Slice only: Select Gantry Sub Type (SLIM GANTRY or NORMAL GANTRY) to switch between LightSpeed HP60 16 slice and BrightSpeed Elite Select. LFC is NOT required.

- 10.) Select the **NETWORK** button to display the Network Settings Screen as shown in [Figure 5-14](#).

Figure 5-14 Network Settings Screen

- 11.) Configure Network Settings

Comment: This screen provides the ability to declare the CT system on a hospital network. Key information such as Host Name, IP Address, Net Mask (for CT systems on a subnet) must be obtained from the hospital network administrator.

See [Chapter 7](#) for more information and complete details of setting the Hospital/System Network Configuration.

- a.) Enter the Suite Name.

The Suite Name must start with a letter, followed by three alphanumeric characters. Total must be four characters long. The name of the OC interface will be <Suite Name>, within the scanner's subnet.

Typically, you should use su01 or ct01 ("su" and "ct" must be lowercase), unless the customer prefers a different suite name.

- b.) Enter the hospital provided Host Name.

The Host Name identifies the network hostname and AE Title of the CT system to the hospital's network. The Host Name:

* **MUST NOT** exceed 16 Characters.

* **MUST** only contain the following characters: a through z, 0 through 9, - and _.

Comment: The Host Name is typically stmary or ct01.

- c.) Enter the hospital provided IP Address for the system.

- d.) Enter the hospital provided Net Mask (if the CT system is on a subnet).

- e.) Enter the Broadcast Address.

Comment:

The Broadcast Address should be the same as the IP Address except for the bits of the host id portion (last digit group) set to 1's or 0's depending on the configuration of the network. The standard default is 1's, but older Sun OS machines used 0's.

For example:

If the IP Address is 192.100.9.17, the Broadcast Address should be 192.100.9.255 if the network is configured to use 1's to specify the Broadcast Address.

If the network contains genesis based scanners or other Sun OS 3.5 or 4.1 computers, the Broadcast Address should be 192.100.9.0.

- f.) Enter the hospital provided Default Gateway IP Address in the Default Gateway field (if applicable). If the site network does not use a default gateway, leave the field blank.
- g.) Select NIS (a.k.a. Yellow Pages database) Advanced Option—only if requested by the hospital network administrator—as follows:
- * Select ADVANCED OPTIONS button on the Network Settings Screen.
 - * Select USE NIS button.
 - * Enter the hospital provided Domain Name (from the hospital network administrator).
 - * Review all pages to be sure the information entered is correct before proceeding to the next step.

- 12.) Select ACCEPT on the System Configuration Screen.

Comment:

The system loads the CT Application Software, OS patches, kernel changes and configures the system on the OC.

The loading process takes approximately 15 minutes. While the load is going on, the results are displayed in a shell window that closes when the loading process is complete. All the window output is logged to a file named: _____

- 13.) When the loading process and configuration changes are complete, the system displays a prompt to reboot. Click on YES. (See [Figure 5-15](#).)



Figure 5-15 Reboot Screen

- 14.) The system will automatically login as ctuser after the reboot. Select OK on the Autostart Disabled popup message.

2.7 Check/Set Date and Time

Tools Required: None

If date and time need to be corrected:

- 1.) Open a Unix Shell window to check the current date.
Type **{ctuser@hostname} date ENTER** at the prompt.

- 2.) If the date is correct, proceed to the next section; if it is incorrect, continue this procedure.
- Note: You must set the date and time on the Host computer with the Application software down.
- 3.) Open a Unix Shell window and login as root:
 - a.) Type `{ctuser@hostname} su -` ENTER at the prompt.
 - b.) Type the password: `#bigguy`
 - 4.) Set the date and time.
 - a.) Type `{root@hostname}# setdate` ENTER
 - b.) Follow the instructions of the individual time-entry prompts, which will appear in the following sequence:
 - * Note: Type "q" at any time to quit. Press ENTER to continue.
 - * Note: TO BE ACCURATE, this tool will prompt you to enter the "Second." Watch your clock or PC carefully to enter the proper value, and hit [ENTER] at the right second to set the accurate time. Enter to proceed. Press ENTER to continue.
 - * Enter the current Year (1980 - 2030) [2007]:
 - * Enter the current Month (1-12) [04]:
 - * Enter the current Day (1-30) [14]:
 - * Enter the current Hour (Military Time) (0-23) [18]:
 - * Enter the current Minute (0-59) [13]:
 - * Enter the current Second (0-59) [00]:
 - * Updating the time on the OC and DARC, Please Wait...
 - * PING darc (172.16.0.2) 56(84) bytes of data.

2.8 Save System State


- 1.) Insert the DVD into the Peripheral Tower DVD drive.
 - 2.) Select: SERVICE DESKTOP.
 - 3.) If reloading software, select: UTILITIES. If upgrading from earlier version software, select: PM.
 - 4.) Select: SYSTEM STATE.
- Note: System State Save may be under Utilities or PM.
- 5.) Select ALL to save all data.
 - 6.) Click SAVE.
 - 7.) If applicable, click OK when the following message appears: System State Media Status: Please insert a DVD or MOD into the drive and press Save again.
 - 8.) When completed select DISMISS.
 - 9.) Label and date the disk including the suite name.
 - 10.) Close the Service Desktop window at the upper left corner of the screen.

2.9 Applications Start-Up

Open the Console shell window, and type: `st` ENTER. The applications desktop appears on the monitor.


2.10 Console Boot-up Flow Chart

Figure 5-16 GOC4/LCGOC Boot-up Flow Chart

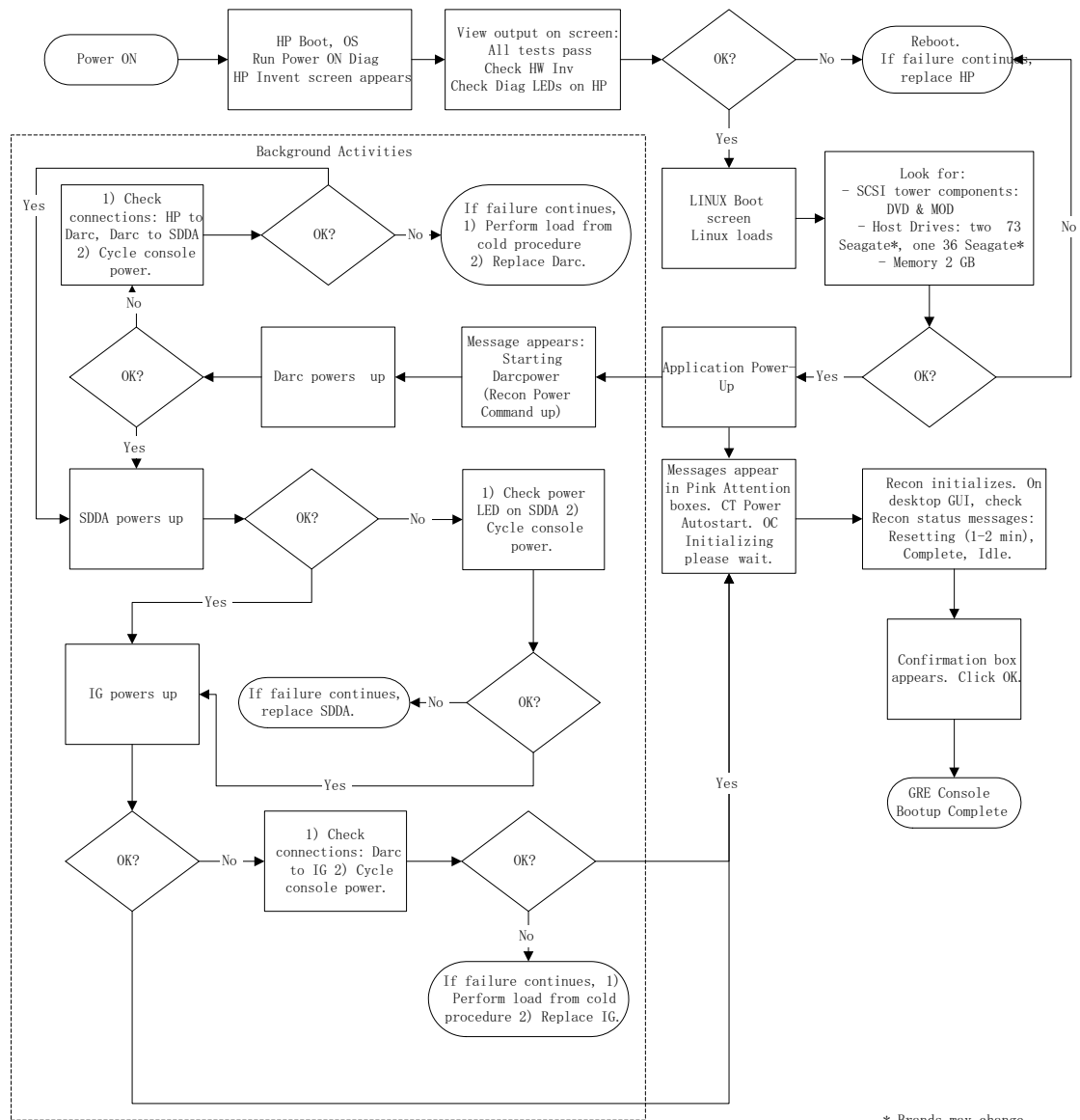

* Brands may change

Figure 5-17 All-In-One Console Boot-up Flow Chart

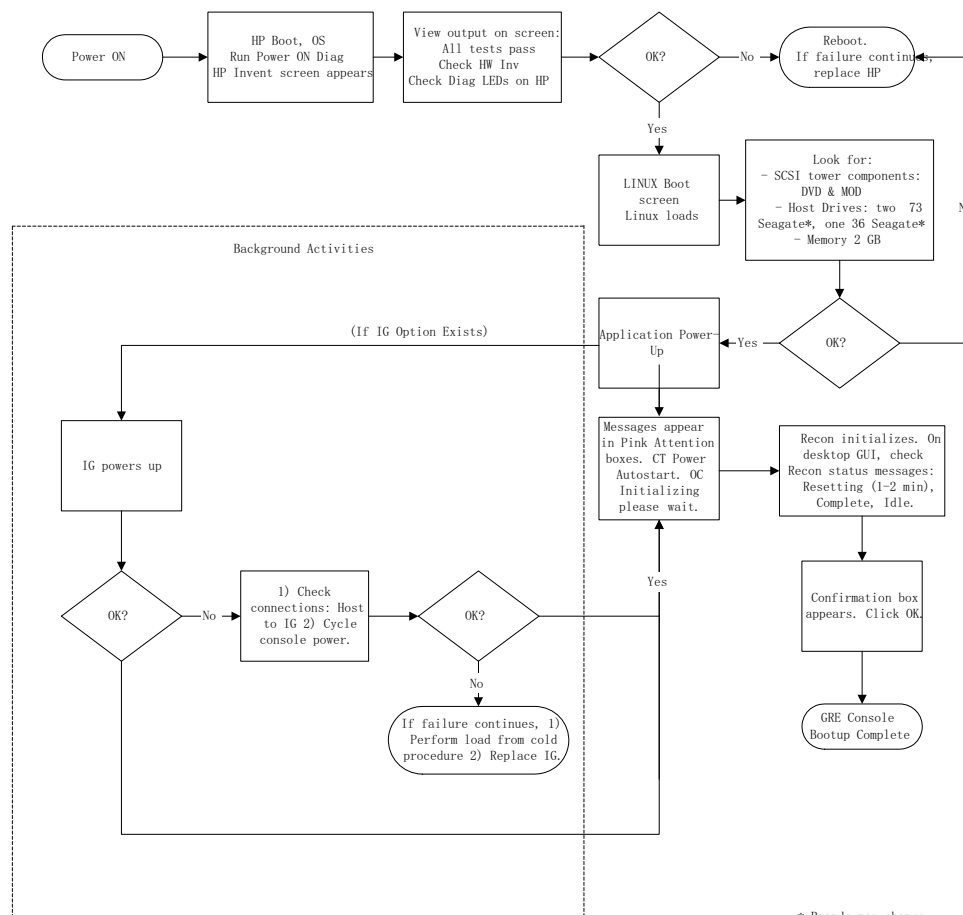

2.11 Adjust Monitor

- 1.) If you are not on the Service Desktop, click on the SERVICE DESKTOP icon. 
- 2.) Click on the IMAGE QUALITY icon. 
- 3.) Select INSTALL SMPTE IMAGE and wait approximately 3-4 minutes for SMPTE image to install. (When complete the following message will display: SMPTE and QA images have been successfully copied)
- 4.) Press ENTER to exit the Service Desktop.
- 5.) Click the IMAGEWORKS icon. 
- 6.) Display the SMPTE pattern. Use the browser to select Exam 1000, which contains the SMPTE pattern, and enlarge the image to full screen display.

-
- 7.) Select Viewer.
 - 8.) Select 1:1 format.
 - 9.) Increase the monitor's contrast to maximum.
- Note: Adjust monitor contrast until the operator sees the anatomical structure (window raster)
- 10.) Increase the Brightness to maximum.
 - 11.) Decrease the Brightness, until the raster just fades into, and matches, the monitor screen background. At this point, the 5% and 95% patches should be just visible.
 - If additional tweaking is required to attempt to match the monitor image to the filmed image, use only the brightness control.
 - If the LCD image exhibits any tearing or smearing of the alphanumeric characters, then reduce the contrast setting slightly until the tearing/smearing is just eliminated. The optimum setting for contrast is the highest setting that does not cause tearing/smearing of the alphanumeric characters.
- You should always finish up by displaying and filming images of anatomy (typical heads and bodies), and asking the technologist to compare the LCD image to the film image.

Section 3.0 Table Gantry Integration

3.1 Introduction

Use these procedures to functionally check every part of the table/gantry subsystem.

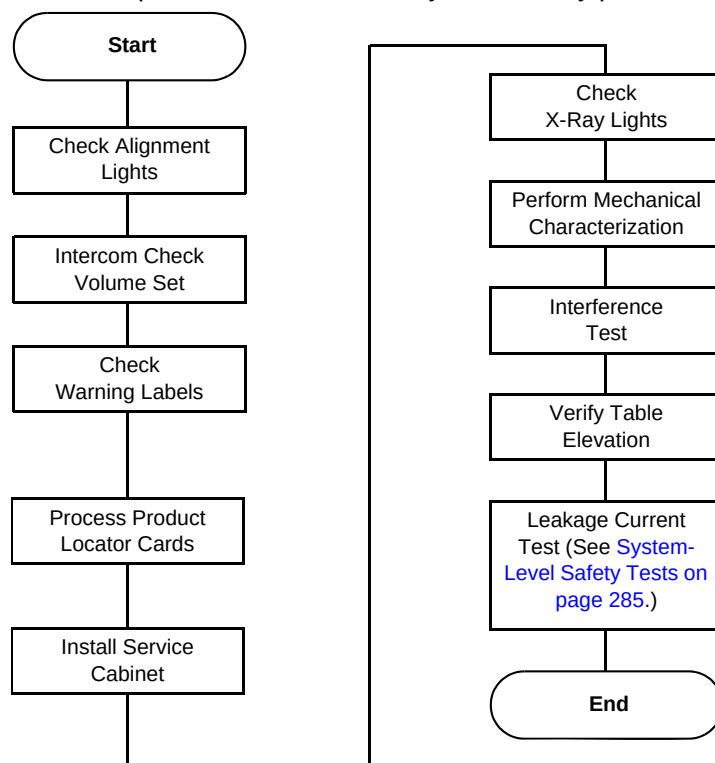


Figure 5-18 Table Gantry Integration Process Overview

Required Tool

- Multimeter

3.2 Check Alignment Lights

3.2.1 Room Light Adjustment

Adjust the scan room lights to normal customer operating levels.

3.2.2 Turning the Alignment Lights ON



CAUTION Verify all personnel have cleared the system. The Gantry rotates during this check.



- 1.) Turn ON the AXIAL DRIVE ENABLE and HVDC ENABLE switches (located on the service switch board).
- 2.) Turn on the alignment light switch on the gantry control panel. The gantry will rotate and the alignment lights will turn ON.



CAUTION LASER EYE INJURY!

NEVER STARE DIRECTLY INTO THE LASER BEAMS WHEN YOU OPERATE THE ALIGNMENT LIGHTS. STARING INTO THE BEAMS CAN CAUSE PERMANENT EYE DAMAGE.

3.2.3 Internal Axial Lights

- 1.) Place a sheet of plain white paper over the output port of each light.
- 2.) Verify that the two laser lines coincide and appear as a single line.

Note: GE designed the internal axial lasers on the current CT system to shine *down* on the collimator. Do NOT adjust the internal alignment lights at this time. The tomographic plane tests use the QA phantom to check the internal axial lasers alignment to the collimator.

3.2.4 External Axial to Internal Axial Distance

Note: Ensure that cradle is level.

- 1.) Raise the table to its highest elevation.
- 2.) Extend the cradle until you see both the internal and external laser lights shining on the cradle.
- 3.) Place a metric rule on the right edge of the cradle, and measure the distance from the internal axial laser line to the external axial line. Verify this distance equals 240.0 mm +1.0 mm.
- 4.) Place the rule on the left edge of the cradle, and measure again.
- 5.) Leave the cradle in its current position, and lower the table to the minimum elevation.
- 6.) Measure the distance between the internal and external lights on both edges of the cradle, as above. Verify the distance remains equal to 240.0 mm +1.0 mm.

3.2.5 Coronal Lights

- 1.) Place a sheet of plain white paper at the left side of the patient opening, in front of the coronal laser light. Verify the two coronal lines coincide.
- 2.) Move the paper to the right side of the patient opening. Verify the two coronal lines coincide.
- 3.) Place the paper in the center of the gantry opening. Use a level to verify that the coronal lines are horizontal.

3.2.6 Turn Lights OFF

Press the alignment light button on the gantry control panel, again, to turn the lights OFF.

3.3 Autovoice/Intercom Check

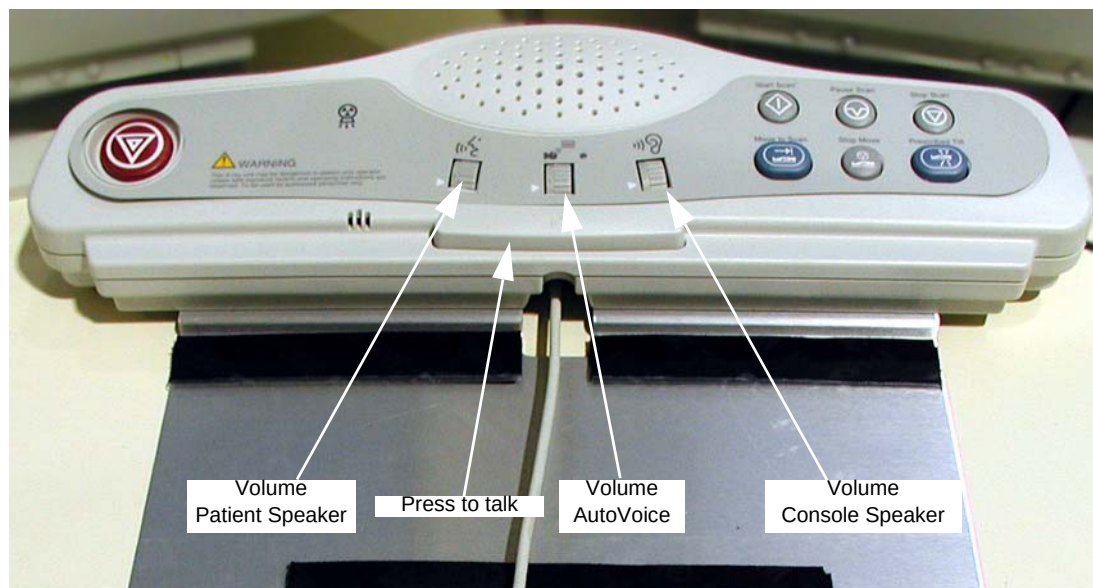


Figure 5-19 SCIM Volume Controls

3.3.1 Requirements

Two people are required to complete this procedure.

3.3.2 Patient Speaker

- 1.) To adjust the volume of the patient speaker in the table, adjust the left-most volume thumb wheel on the SCIM while speaking into the console microphone. (Press the bar on the SCIM to talk; release the bar to listen.)
- 2.) The patient should be able to clearly hear the operator.

3.3.3 Operator Console Speaker

To adjust the console speaker volume:

- 1.) Have an assistant speak into the gantry microphone.
- 2.) Adjust the SCIM console volume knob until you can clearly hear the patient.

3.3.4 Autovoice Volume

- 1.) On the Scan Desktop, select PROTOCOL MANAGEMENT.
- 2.) Select AUTO VOICE RECORD.
- 3.) Click the 3.4 button, to the right of "FF2. Inspiration".
- 4.) Click the PLAY button, to play the Inspiration AutoVoice message.
- 5.) Adjust the center volume thumb wheel while Autovoice is playing, to set the volume for the gantry speaker.
- 6.) Repeat steps 4 and 5 as necessary to achieve satisfactory volume.
- 7.) Select DONE, then select QUIT.

Note: If a satisfactory volume can not be achieved, refer to the system service manual and review the intercom module setup procedure.

3.4 CT System X-Ray ON Indicators, Cautions & Warning Labels

3.4.1 Check And Install System Warning Labels

All labels are installed in English and present on PDU, Console, Table, Gantry and Accessories. Add the labels listed below ([Table 5-2](#) and [Section 3.5](#)) for the appropriate language for the country in which this system is installed. Additionally, apply any other warning labels if present, on equipment where appropriate.



NOTICE Do not cover English labels already on the system.

Subsystem	Component	Label(s)
Console	SCIM	SCIM overlay warning label
	Keyboard	Function key overlay label
Gantry	Scan Window	Laser warning label
	Laser Window	Laser warning label
	Front Cover	Laser warning label Information labels
	System GIB	System Global Installation Base (rating) label
Table	Top Side Cover	Pinch Hazard warning label - each side of cover
NGPDU	Front Cover	Emergency OFF label Gantry Enable label Power ON label
Accessories	Table Foot Extender	Warning label
	Coronal Head Holder	Warning label
	Sagittal Head Holder	Warning label
	Accessory Tray	Warning label
	IV Pole	Caution label

Table 5-2 System Warning Labels

3.4.2 Documentation - Verification

When finished update GE form e4879 (US Only) and the installation completion form that all appropriate language labels were installed and present.

3.5 Check Warning Labels

Note: The labels on the system and the system manuals must comply with the country law, as listed in Direction 5221102-1EN (found in the keyboard collector kit shipped with the system) regardless of the user interface (UI) language that is chosen. Compliance to the law must be completed prior to releasing the system to the customer.

Note: Do not cover English labels already on the system.

3.5.1 On SCIM

- 1.) Make sure the SCIM overlay warning label appears in the correct location on the SCIM.
- 2.) Record this information on Form 4879. For more information about this form, see [Section 8.0](#) of Chapter 4.

3.5.2 On Gantry

- 1.) Check that all laser warning labels are present on the gantry near the laser opening.
- 2.) There should also be warning labels on the lower right side of the gantry front cover.
- 3.) Record this information on Form 4879. For more information about this form, see [Section 8.0](#) of Chapter 4.

3.5.3 On Laser

- 1.) Make sure all laser warning labels appear in the correct location on the outside of the gantry.
- 2.) Obtain and install replacements for any missing labels.



Figure 5-20 Laser Warnings and Precautions

3.6 Process Product Locator Cards

- 1.) Collect the product locator cards shipped with the system. There should be approximately 28 product locator cards with the average system.
- 2.) Update the online product locator web site with the required hospital information.
- 3.) Confirm that the serial numbers on the cards shipped with the system match those found on the web site for that GON number.
- 4.) Update as required. Place the cards in a plastic bag, then place them in the service cabinet.

3.7 Verify Service Cabinet Installation

Verify that the service cabinet is installed and that the shipped service materials are in the cabinet. If the cabinet was not installed, install it at this time by following the Service Cabinet Installation procedure in [Section 5.6](#) of Chapter 2.

3.8 Check X-Ray Lights

Perform several scans following the steps below. Verify that the X-ray ON lights are ON during the scans. When done, check the boxes in [Table 5-3](#).

- 1.) Make sure the axial drive enable and HVDC enable switches are ON.
- 2.) If you are not on the Service Desktop, click on the Service Desktop icon.
- 3.) Select DIAGNOSTICS.
- 4.) Select DIAGNOSTIC DATA COLLECTION.
- 5.) Set the scan time to 2.00.
- 6.) Set the kV to 80.
- 7.) Set the mA to 40.
- 8.) Press ACCEPT RX.
- 9.) Press START SCAN button when flashing.
- 10.) Record the above information on Form 4879. For more information about this form, see [Section 8.0](#) of Chapter 4.

Light On?	Warning Light Locations
☐	SCIM
☐	Gantry Front
☐	Gantry Back (Use a mirror to view)
☐	Room Light (outside of the room)

Table 5-3 X-ray Light Chart

3.9 Mechanical Characterization

The relationship of table height to ISO center and internal-to-external landmarks must be characterized for proper interference matrix functionality.

Note: Do NOT perform tilt characterization.

3.9.1 Alignment Light Characterization

- 1.) Start the Mechanical Characterization tool from the Calibration tab on the Common Service Desktop.
- 2.) Select the CHARACTERIZE ALIGNMENT LIGHTS button from the interface.
- 3.) Follow the on-screen instructions.

3.9.2 Table Height Characterization

- 1.) Select the CHARACTERIZE TABLE HEIGHT button from the interface.
- 2.) Follow the on-screen instructions.

Note: If the table height is less than 21mm or greater than 25mm, relative to ISO, you must adjust the table height using the table leveling pad and adjusters. Raise or lower all four adjusters equally to achieve desired results. Note down the value of distance "V" which will be used in Section 3.10.

3.10 Interference Test

PREREQUISITES

- Be sure that the System State was restored from **DVD** per Section 2.0.
- Reset the hardware to download the new characterization values before performing the table/gantry interference tests in this section.

CONFIGURATION

- 1.) Table flashed with latest SW
- 2.) Table must have elevation and cradle characterized
- 3.) Table must be mechanically aligned to gantry
- 4.) Table must have had the table-gantry characterization completed
- 5.) Verify the Table extender is installed.

TEST OUTLINE

The following tests verify the proper tilt and table interference matrix on the H-Power gantry.

- [Verify Table Elevation – Section 3.10.1](#)
- [Position Tilt, Move Table to Interference Limit – Section 3.10.2.1](#)
- [Position Table, Move Tilt to Interference Limit – Section 3.10.2.2](#)
- [Tilt Limits When Table Below Scan Plane Lower Limit – Section 3.10.2.3](#)

REQUIREMENTS

The following requirements are tested in this series of tests.

- 1.) No motion shall cause the table to hit the gantry (or gantry to hit the table)
- 2.) Requirement #1 shall include the use of the table extender.
- 3.) No tilt motion shall cause the gantry tilting frame to touch the stationary base covers for any tilt angle.

LIMITATIONS

These requirements will only be met when the table is NOT in service mode.

INTERPRETING TEST RESULTS

If test results indicate that elevation and/or tilt display readings fail to meet specifications, DO NOT adjust the limit switches. Instead, re-characterize and/or adjust tilt speed. Refer to Appendix C. Perform elevation and cradle first, then repeat the test. If it still fails, perform tilt.

3.10.1 Verify Table Elevation

Note: "V" means distance from table height to ISO.

#	TEST	EXPECTED RESULTS
3.10.1-1	Move the cradle to home position. Push the table down gantry push-button to lower the table to the minimum height.	Elevation Display should read $550.0 + V \pm 3 \text{ mm}$

Table 5-4 Table Elevation Tests

#	TEST	EXPECTED RESULTS
3.10.1-2	Push the table up gantry push button to raise the table to the maximum height. If the mechanical alignment of the table/gantry is not correct, as is often the case during manufacturing staging, this value may be out of range. Most of the following tests will still be valid: Those that might have some variation are indicated in bold type.	Elevation Display should read $V \pm 3$ mm. If mechanical alignment of the table/gantry is not correct, this value can be as low as 0 mm and as high as 40 mm.

Table 5-4 Table Elevation Tests

3.10.2 Tests

3.10.2.1 Position Tilt, Move Table to Interference Limit

The following tests verify the table interference limits at different tilt locations.

- Note:
- "I" means top of gantry tilts toward the table base
 - "S" means top of gantry tilts away from the table base.
 - "V" means distance from table height to ISO.

Important: For all tests, make sure there is 2.5 cm of clearance between the gantry and table. Also make sure that the Gantry Tilting frame covers do not touch the stationary base covers, for all tilt angles used in this test.

#	TEST	EXPECTED RESULTS
3.10.2.1-1	Move cradle to home position and set internal landmark.	Cradle position on display should read 0.0
3.10.2.1-2	Raise table height to maximum height, set the internal landmark, move cradle into gantry 1m.	Table elevation on display should read $V \pm 3$ mm. (This value will be the same as in 3.10.1-2 above.) Cradle position on display should read 1000.0 mm.
3.10.2.1-3	Tilt the gantry to I30.0, then lower table until motion stops.	Table elevation on display should read 100 ± 3 mm. Tilt display should read I30.
3.10.2.1-4	Tilt the gantry to I23.0, then lower table until motion stops.	Table elevation on display should read 155 ± 3 mm. Tilt display should read I23.
3.10.2.1-5	Tilt the gantry to I20.0, then lower table until motion stops.	Table elevation on display should read 170 ± 3 mm. Tilt display should read I20.
3.10.2.1-6	Raise the table elevation to maximum height.	Table elevation on display should read $V \pm 3$ mm. (This value will be the same as in 3.10.1-2 above.)
3.10.2.1-7	Tilt the gantry to S30.0, then lower table until motion stops.	Table elevation on display should read 156 ± 3 mm. Tilt display should read S30.
3.10.2.1-8	Tilt the gantry to S23.0, then lower table until motion stops.	Table elevation on display should read 182 ± 3 mm. Tilt display should read S23.
3.10.2.1-9	Tilt the gantry to S20.0, then lower table until motion stops.	Table elevation on display should read 192 ± 3 mm. Tilt display should read S20.
3.10.2.1-10	Raise the table to 154 mm.	Table elevation on display should read 154mm.

Table 5-5 Position Tilt, Move Table to Interference Limit Tests

#	TEST	EXPECTED RESULTS
3.10.2.1-11	Tilt gantry to S30 and verify the table height can be adjusted from 154 to 25 mm. (This value will be the same as in 3.10.1-2 above.)	Tilt display should read S30. Table lower limit should be 154 ± 3 mm. Upper table limit should be $V \pm 3$ mm. (This value will be the same as in 3.10.1-2 above.)
3.10.2.1-12	Set the table height to 98 mm.	Table elevation on display should read 98mm.
3.10.2.1-13	Tilt gantry to I30 and verify the table height can be adjusted from 98 to 25 mm. (This value will be the same as in 3.10.1-2 above.)	Tilt display should read I30. Table lower limit should be 98 ± 3 mm. Upper table limit should be $V \pm 3$ mm. (This value will be the same as in 3.10.1-2 above.)

Table 5-5 Position Tilt, Move Table to Interference Limit Tests

3.10.2.2 Position Table, Move Tilt to Interference Limit

The following tests verify the tilt interference limits at different table heights.

- Note:
- "I" means top of gantry tilts toward the table base
 - "S" means top of gantry tilts away from the table base.
 - "V" means distance from table height to ISO.

Important: For all tests, make sure there is 2.5 cm of clearance between the gantry and table.

#	TEST	EXPECTED RESULTS
3.10.2.2-1	Move cradle to home position and set internal landmark. Set gantry tilt to zero.	Cradle position on display should read 0.0. Gantry tilt on display should read 0.0.
3.10.2.2-2	Raise table height to maximum height, set the internal landmark, move cradle into gantry 1m.	Table elevation on display should read $V \pm 3$ mm. (This value will be the same as in 3.10.1-2 above.) Cradle position on display should read 1000.0 mm.
3.10.2.2-3	Lower table until height is 115mm. Tilt the gantry top away from the table ("S") until it stops.	Table elevation on display should read 115 mm. Tilt display should read $S30 \pm 0.5^\circ$.
3.10.2.2-4	Tilt the gantry top toward the table ("I") until it stops.	Table elevation on display should read 115 mm. Tilt display should read $I28 \pm 0.5^\circ$.
3.10.2.2-5	Tilt the gantry to 0. Lower table until height is 200mm. Tilt the gantry top away from the table ("S") until it stops.	Table elevation on display should read 200 mm. Tilt display should read $S17 \pm 0.5^\circ$.
3.10.2.2-6	Tilt the gantry top toward the table ("I") until it stops.	Table elevation on display should read 200 mm. Tilt display should read $I14.5 \pm 0.5^\circ$.
3.10.2.2-7	Tilt the gantry to 0. Lower table until height is 210 mm. Tilt the gantry top away from the table ("S") until it stops.	Table elevation on display should read 210 mm. Tilt display should read $S13.5 \pm 0.5^\circ$.
3.10.2.2-8	Tilt the gantry top toward the table ("I") until it stops.	Table elevation on display should read 210 mm. Tilt display should read $I20 \pm 0.5^\circ$.

Table 5-6 Position Table, Move Tilt to Interference Limit Tests

3.10.2.3 Tilt Limits When Table Below Scan Plane Lower Limit

The following tests verify the table and tilt interference limits when the table height is below the scan plane.

- Note:
- "I" means top of gantry tilts toward the table base
 - "S" means top of gantry tilts away from the table base.
 - "V" means distance from table height to ISO.

Important: For all tests, make sure there is 2.5 cm of clearance between the gantry and table.

#	TEST	EXPECTED RESULTS
3.10.2.3-1	Set gantry tilt to zero. Move cradle to home position, lower the table all the way, and set the internal landmark.	Cradle position on display should read 0.0 Gantry tilt on display should read 0.0. Table height should read $550.0 + V \pm 3$ mm. (This value will be the offset from 575 mm by the same amount as the upper table limit varies from 25 mm in 3.10.1-2 above.)
3.10.2.3-2	Tilt the gantry forward and backwards and verify the following tilt limits: S4.0 and I30.	Gantry tilt on display should read $S4.0 \pm 0.5^\circ$. Gantry tilt on display should read $I30.0 \pm 0.5^\circ$.
3.10.2.3-3	Set Gantry tilt to 0 and set the internal landmark. Then, using the gantry push-buttons, move the cradle in towards the gantry until it is stopped.	Gantry tilt on display should read 0.0. Cradle should stop at $91 \text{ mm} \pm 6 \text{ mm}$ from the home position. (This value will not be valid if the table/gantry mechanical characterization is not correct as described in 3.10.1-2 above. Insure that the front of the table is at least 25 mm from the gantry.)
3.10.2.3-4	Set gantry tilt to zero. Move cradle to home position.	Cradle position on display should read 0.0 Gantry tilt on display should read 0.0.
3.10.2.3-5	Raise the table to a height of 386 mm and verify the following tilt limits: S10 and I30.0.	Table height should read 386 mm. Gantry tilt on display should read $S10.0 \pm 0.5^\circ$. Gantry tilt on display should read $I30.0 \pm 0.5^\circ$.
3.10.2.3-6	Set Gantry tilt to 0. Set the internal landmark. Then, using the gantry push-buttons, move the cradle in towards the gantry until it is stopped.	Gantry tilt on display should read 0.0. Cradle should stop at $142 \text{ mm} \pm 10 \text{ mm}$ from the home position.
3.10.2.3-7	Set gantry tilt to zero. Move cradle to home position.	Cradle position on display should read 0.0. Gantry tilt on display should read 0.0.
3.10.2.3-8	Raise the table to a height of 242 mm and verify the following tilt limits: S18.5 and I30.0.	Table height should read 242 mm. Gantry tilt on display should read $S18.5 \pm 0.5^\circ$. Gantry tilt on display should read $I30.0 \pm 0.5^\circ$.
3.10.2.3-9	Set Gantry tilt to 0. Raise the table to 210 mm. Then, using the gantry push-buttons, move the cradle in towards the gantry until it is stopped.	Gantry tilt on display should read 0.0. Cradle should go all the way through the gantry bore to the full-extended position.
3.10.2.3-10	Set gantry tilt to zero. Move cradle to home position, lower the table all the way and set the internal landmark.	Cradle position on display should read 0.0 Gantry tilt on display should read 0.0. Table height should read $550.0 + V \pm 3$ mm. (This value will be the offset from 575 mm by the same amount as the upper table limit varies from 25 mm in 3.10.1-2 above.)
3.10.2.3-11	Tilt the gantry top toward the table to a tilt of 30 degrees	Display should read I30.

Table 5-7 Tilt Limits When Table Below Scan Plane Lower Limit Tests

#	TEST	EXPECTED RESULTS
3.10.2.3-12	With the table down all the way, move the cradle in until it stops.	Cradle position should be 430 ± 10 mm.
3.10.2.3-13	Move the table to the home position, raise the table to a height of 386 mm, set the internal landmark, and then move the cradle in until it stops.	Cradle position should be 413 ± 10 mm.
3.10.2.3-14	Move the table to the home position, raise the table to a height of 100 mm, then move the cradle in.	Cradle should go all the way through the gantry bore to the full-extended position.

Table 5-7 Tilt Limits When Table Below Scan Plane Lower Limit Tests

Chapter 6

Image Quality



NOTICE Potential for Data Loss and/or Equipment Damage

To prevent potential data loss, please do the following:

- When instructed, record data collected from the procedures in this chapter on GE Form e4879. For more information about this form, see [Section 8.0](#) of Chapter 4.
- Only use the Installation manual that arrives with your system for installation. Any other revisions of this manual may not exactly match your system.

Section 1.0

Introduction

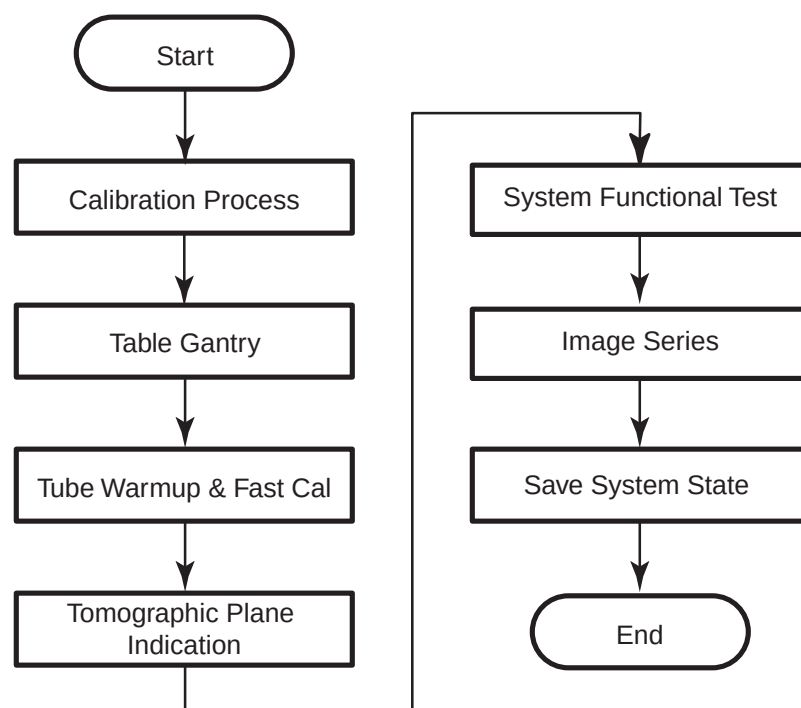


Figure 6-1 System Test Overview

Section 2.0

Calibration Process

2.1 Reference Procedures

Do not perform these procedures ([Scanning with Service Protocols](#) or [Center Phantom](#)) until instructed to do so in other sections of this chapter.

2.1.1 Scanning with Service Protocols

Locate the Manufacturing and Installation protocols under Infant area 20.

Note: **Manufacturing and Service share this Protocol list. Different product option offerings also use this list. Carefully follow the scan section instructions, and verify you acquired the images with the correct technique before filling out the data sheet. Otherwise you may troubleshoot an image problem that only exists because you used the wrong technique.**

- 1.) Select the NEW PATIENT icon on the left monitor
- 2.) Enter a Patient ID (e.g., **getest**).
- 3.) Click on the box labeled PEDIATRIC.
- 4.) Select a service protocol from the list, to display the corresponding view edit screen.
Optional Method: Enter the protocol number into the Protocol Number Field on the Exam Rx Screen.

2.1.2 Center Phantom

- 1.) Select SCANNER UTILITIES icon on the left monitor
- 2.) Select CENTER PHANTOM.
- 3.) Follow the on-screen procedures.
- 4.) The phantom center spec is $\pm 5\text{mm}$.
- 5.) Select QUIT, when the phantom is within specification.
- 6.) Level the phantom both front to back and side to side. (use a 6" level)

2.2 Prepare the QA Phantom

- 1.) Locate the multi-language sticker packet in the QA phantom shipping box.
- 2.) Attach the sticker with the customer's language to the face of the phantom hanger bracket.

2.3 Calibration Process Introduction

If your system has a factory supplied state DVD-RAM, you used it to load the system calibration files during the Restore System State ([2.3 on page 199](#) of this manual).

Section 3.0

Table/Gantry Alignment Procedure

3.1 Time & Personnel

Required Persons	Preliminary Reqs	Procedure	Finalization
2 (FE or mechanical supplier)		45 minutes labor on-site	

3.2 Tools and Test Equipment

- 1 mm wire

3.3 Preparation

- All table mechanical alignment procedures completed.
- The table perpendicular alignment test passed.
- Table anchors are in place and within specification.
- The table is level in all directions.

3.4 Procedure

TABLE GANTRY PREP

- 1.) Check that the table cradle is level in all directions. Correct, if necessary.
- 2.) Drive the table to its highest elevation ISO with the phantom holder removed.
- 3.) Check the scan window for proper installation.

VERIFY TABLE DRIVE CONSISTENCY

- 4.) Drive the table cradle in and out five times to seat the rollers.

CRADLE SETUP

- 5.) Turn on the alignment lights.
- 6.) Advance the end of the cradle to the black dot on cradle.
- 7.) Tape a 100 mm section of 1 mm wire on the cradle that aligns with the white cradle center line.
- 8.) Using the gantry keypad, set an internal landmark, and then advance the cradle 1000 mm.
- 9.) Tape second 100 mm section of 1 mm wire on the cradle that aligns with the alignment light.

SCANNING SETUP

- 10.) From the application screen Select *NEW PATIENT*.
 - a.) Fill out patient ID: GE Test
 - b.) Name: Alignment
- 11.) From the Protocol screen:
 - a.) Select *SERVICE*,

- b.) Select *IMAGE QUALITY*,
- c.) Select *PERPENDICULAR ALIGNMENT*.
- 12.) The red boxes should disappear from the screen. If not, reset an internal landmark.
- 13.) Select *CONFIRM*, then press the START SCAN button when lighted.

IMAGE REVIEW

- 14.) On the Service screen, select *IMAGE WORKS*
 - a.) Locate the scanned examination in the Examinations column.
 - b.) Highlight the Alignment scans
 - c.) Select *VIEWER*.
 - d.) Select *FORMAT*, and select the two-in-one format horizontal display view.
- 15.) Click on image 1 and select the grid. With the grid and image displayed, visually compare image 1 to image 2.
 - For close inspection, you may need to use the zoom function to see a difference.
 - Visually compare image 1 and image 2 to verify the centering wire appears in the center of the grid. As shown on the screen, the wire is 1 mm. Use the measure tool to determine the alignment difference. Move the table until both are within ± 2 mm of center.

MOVING THE TABLE

- 16.) The adjustment is likely to require a very small movement. Use a suitable tool to move the table the required distance.
- 17.) Rescan to confirm each movement trial. This can take a few trials to move the table to a position that is within the specification.
- 18.) The plastic accessory edges of the cradle are installed to allow cradle accessories to be used. These edges, although visible, should not be used to determine cradle center. Edge-to-edge difference can be greater than the alignment specification.
- 19.) Repeat above steps until both images are visually aligned on the screen.

FINALIZATION

- 20.) Use a calibrated torque wrench to tighten the anchors to 50 ± 6 N-m (37 ± 5 ft.-lb.). Confirm that the torqued anchor still meets the anchor installation specifications:
 - a.) Maintain 1 full thread of adjustor showing above the lock ring or table base plate.
 - b.) Have not more than 1 in. of anchor showing above the nut. Do not cut off any access.
 - c.) Using a permanent marker, draw a line on the nut and base. Use this line to determine whether the anchor loosened over time.
- 21.) Reinstall all table components removed to access the anchors.

Section 4.0 Tube Warm Up and Fast Cal



- 1.) Select  to warm up the tube.
- 2.) Select FAST CALIBRATION from the Daily Prep menu.

- Note: **Use the default Fast Cal selections determined by the system configuration. (The system defaults to all four kV stations, but you can choose kV stations to calibrate during reconfig.)**
- 3.) Run the selected air calibrations.
 - 4.) When the calibration process completes, click on QUIT.

Section 5.0

Tomographic Plane Indication

- 1.) Place the QA phantom on the phantom holder.
- 2.) Center the Phantom (refer to procedure [2.1.2 on page 228](#))
- 3.) Turn ON the internal alignment lights, and drive the phantom into the gantry opening, until the line on the phantom lines up with the internal laser lights.
- 4.) Verify that BOTH internal axial lasers line up along the line on the QA phantom. If not, check table/gantry, cradle, and/or laser alignment.
- 5.) Center the phantom in the scan plane with the calibration program. (See [2.1.2 on page 228](#) for details on the phantom centering procedure.)
- 6.) Select the service protocol TOMO PLANE INDICATION from MANUFACTURING. (See [2.1.1 on page 228](#) for details on scanning with service protocols.)

or

Manually select the scan parameters in [Table 6-1](#).

Scan Type	kV	mA	SFOV	Thickness	Scan Time	Start Loc.	End Loc.	Algorithm	Interval
Helical	120	160	Small	1.25HQ	1.0sec	I3.0	S3.0	Bone	0.2

Table 6-1 Tomographic Plane Indication Scan Parameters

- 7.) Display the image series, and locate the scan plane indicator, the longest bar in the bar pattern on the right side of the phantom. The right side of the phantom corresponds to the side of the image labeled **L** on the display screen.
- 8.) On the HHS Data Sheet, record the scan location (shown on the image annotation) of the image with the darkest scan plane indicator (darkest long bar).
- 9.) If your system meets all the installation and alignment specifications, the image at scan location zero (S0.0) should contain the scan plane indicator. If scan location S1.0 or scan location I1.0 has the darkest bar, the system still meets the specification. The scan plane deviation should equal $S0.0 \pm 1.0\text{mm}$. If necessary, adjust the internal alignment light position to meet the $S0.0 \pm 1.0\text{mm}$ requirement.
- 10.) Repeat the Tomographic Plane Indication test with the external alignment lights.
 - a.) Use the external alignment light, and press the external landmark.
 - b.) Verify the external light lines up along the black line on BOTH the left and right sides of the QA phantom.
 - c.) The scan plane indication must fall within the $S0.0 \pm 1.0\text{mm}$ specification.
- 11.) Initial below.

Section 6.0 Image Series

IMPORTANT: Run ALL Image Series Tests in Auto Mode. Manual procedures are provided in this section as REFERENCE ONLY.

6.1 Scan Protocol

The person who acquires the image series has the responsibility to review the images and verify they meet the specifications listed on data sheets. Responsibilities also include means and standard deviation measurements and keeping a record of failures that occur during the image series. Unless otherwise stated, use the following scan parameters during the image series acquisition:

- Scan FOV equal to display FOV (Field of View)
- 512x512 matrix size

Note: **Consider any image series scan that does not meet specifications as failing.**

For means and standard deviations, 90% of the slices must pass. Any failure on a particular technique requires at least ten additional slices to evaluate effectively.

Systems with metal-free cradles have a phantom holder with a perpendicular adjustment (Z-axis) knob on it. **Each time you change phantoms**, make sure you use a bubble level, and the Z-axis knob on the phantom holder, to level the phantom.

6.2 Data Recording: Means and Standard Deviation

Any failure on a particular technique requires at least a ten additional slices to evaluate effectively. For means and standard deviations, 90% of the slices must pass.

- Record means to two decimal places, and round to the nearest one-tenth, (one decimal place) when you compare the resulting values to the spec.
- Record standard deviations to two decimal places, then round off to one decimal place, to compare it to the spec.
- Average standard deviations: Use two decimal places to average the values, then round off to one place.

Before you record the means and standard deviations, check the image data sheets to determine whether to average the means and standard deviations, or record them slice by slice. Make sure you record all the required image data on the HHS data sheets.

6.3 Term Definitions and Screens

Xc - Mean CT number for the specified center coordinates of the phantom image.

AvXc - Average Mean CT number for the center of the phantom image: Average the mean CT value for all specified center coordinates of all slices in an exam.

Xo - Mean CT number for the outside of the phantom image: Average the mean CT value for all specified outside coordinates of one slice.

AvXo - Average outside mean CT number for the number of slices in an exam.

AvSDc - Average image noise about the center image coordinate (measured as the standard deviation) of all slices in an exam.

AvSDo - Average image noise (standard deviation) for the outside of a phantom: Average of all outside coordinates of all the slices in an exam.

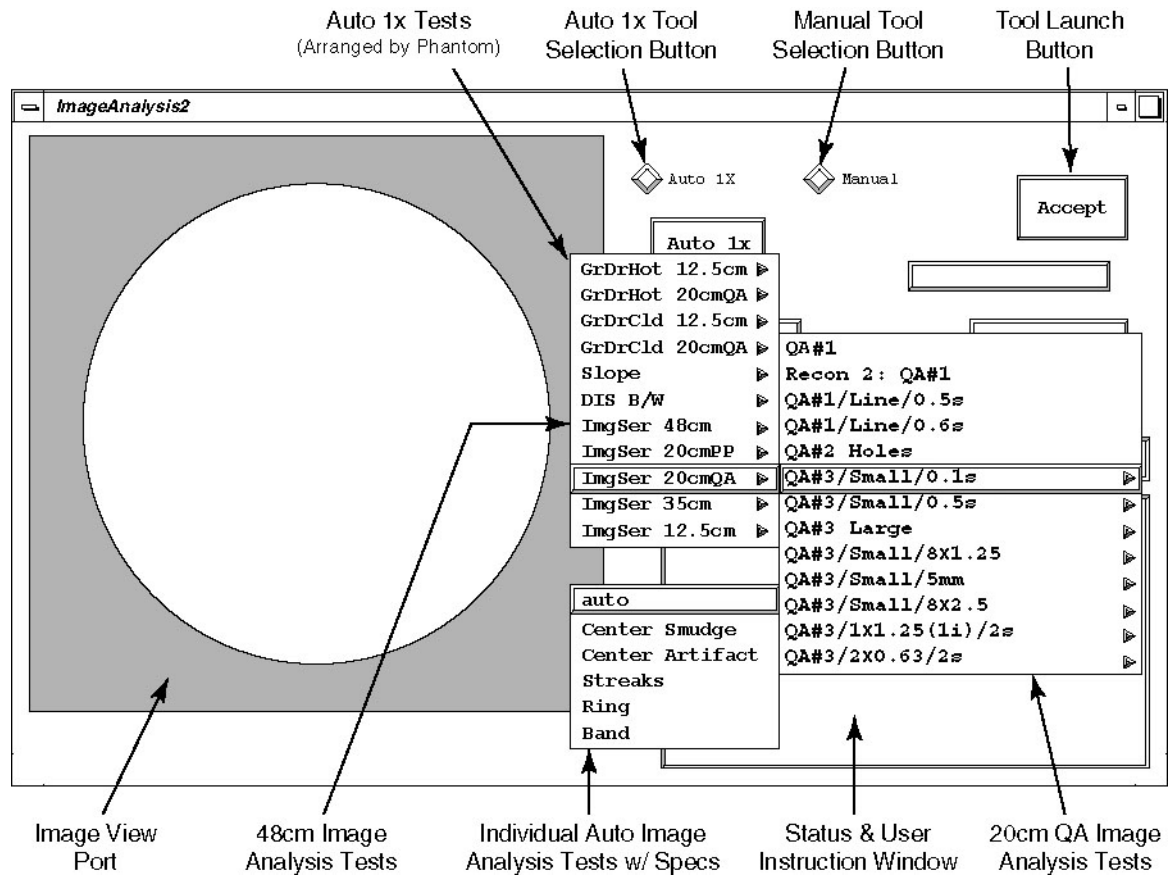


Figure 6-2 Image Analysis Tool User Interface - Auto 1x Test Pull Down Menu

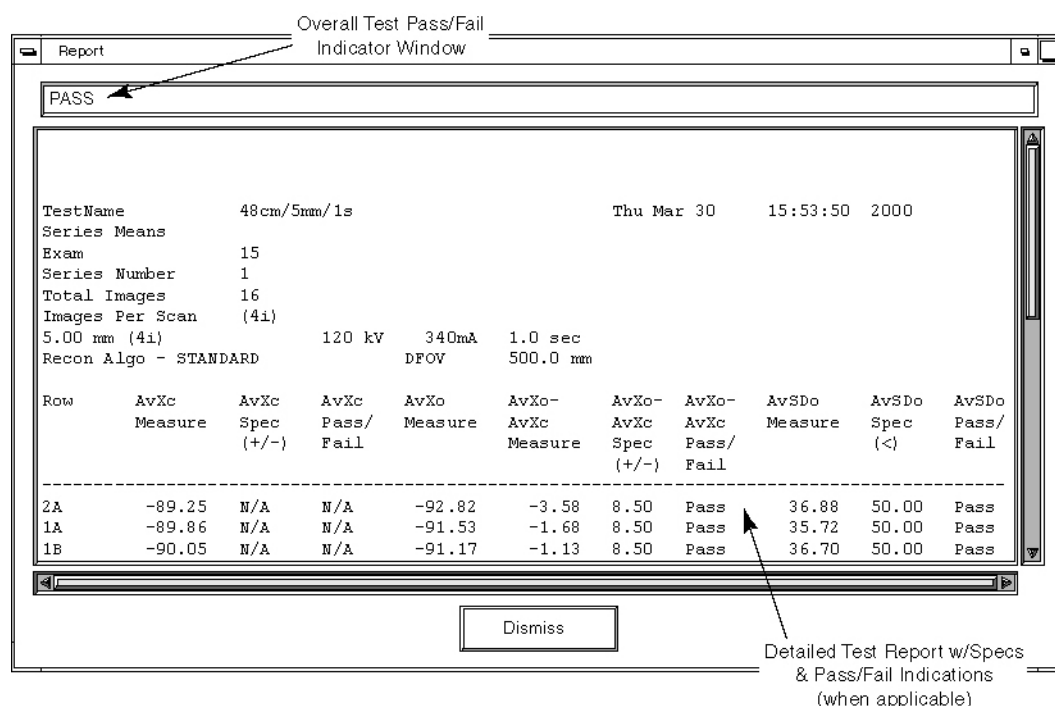


Figure 6-3 Image Analysis Tool User Interface - Test Results Report Window

6.4 48cm Phantom Image Series Image Performance Verification

6.4.1 Acquiring the 48cm Phantom Image Series

- 1.) Mount the Phantom Holder on the head-end of the table.
- 2.) Mount the 48cm Phantom on the Phantom Holder.
- 3.) Align, level, & center the 48cm Phantom.
 - Align phantom using the internal laser lights.
 - Level phantom using bubble level and the Z Axis knob on the Phantom Holder.
 - Center phantom using the CENTER PHANTOM procedure in the left head SCANNER UTILITIES selection and the X and Y Axis knobs on the Phantom Holder.
- 4.) Set up the system to scan a single, 64 image, 48cm Phantom series.

MANUAL SCAN PROTOCOL SETUP

Refer to [Table6-2](#) to manually set-up an Axial scan with the parameters shown.

Note: These parameters may differ from those found using the [Auto Scan Protocol Setup](#) below.

System	Series Description	Scan Type	Start Loc.	End Loc.	Total # of Images	Thick Speed	Interval (mm)	Gantry Tilt	SFOV	kV	mA	Total Exposure Time	DFOV (cm)	Recon Type
16 Slice	48 16x1.25/120kV /200mA/2s	Axial Full 2.0 sec	I9.375	S9.375	64	1.2516i	0.00	S0.0	Large	120	200	8.0 sec.	50.0	Std
4/8 Slice	48 Lg Series 1	Axial Full 1.0 sec.	I7.50	S7.50	16	5.0 4i	0.00	S0.0	Large	120	260	4.0 sec.	50.0	Std

Table 6-2 48cm Phantom Image Series Scan Parameters

AUTO SCAN PROTOCOL SETUP

- a.) On the Exam Rx desktop, select NEW PATIENT.
 - b.) Type the following entries in the two listed Patient Information fields:
 - * Patient ID: **Service**
 - * Name: **48cm Phantom Image Series**
 - c.) From the Protocol Selection display, click on the Service Tab under Anatomical Selector.
 - d.) On the Service Protocol Selection window, click on MANUFACTURING.
 - e.) On the Service Category Protocol List, click on IMAGE SERIES 48CM.
 - f.) **16 Slice**
On the ExamRx protocol parameter display, skip to the 3rd Series.
(Series Description: 48 16x1.25/120kV/200mA/2s)
4/8 Slice:
On the ExamRx protocol parameter display, skip to the 2nd Series.
 - g.) Set internal Landmark.
- 5.) Acquire a single, 64 image, 4 scan, image series of the 48 cm Phantom.

6.4.2 Image Performance Verification

6.4.2.1 Series Means (AvXo-AvXc and AvSDo) & Center Smudge

- 1.) Launch the Image Analysis Tool.
 - a.) From the Global Control Palette, click on the Service Desktop.
 - b.) From the Service Desktop upper navigational bar, click on the IMAGE QUALITY icon.
 - c.) In the left navigational area directory, click on IMAGE ANALYSIS.
- Note: During the Image Analysis process, you will need to switch between the Service Browser screen (to select the image series to analyze) and the Image Analysis screen (to run a particular Auto 1x analysis test). Use ALT/F1 to bring a window to the front. Use ALT/F3 to send a window to the back.
- 2.) Using the AUTO 1X Tool on the Image Analysis screen, generate a Series Means and Center Smudge report for the 1st series of the 48cm Phantom acquisition.
 - a.) From the Service Browser window, select the 48cm Phantom series acquired in Section 6.4.1.
 - b.) On the Image Analysis Tool window, click on the AUTO 1X diamond shaped button. (See [Figure 6-2, on page 234](#).)
 - c.) Click on the AUTO 1X rectangular button, select the IMGSER 48CM test from the pull-down menu, slide cursor to the right and select 16X1.25/120KV/200MA/2S (**16 Slice**) or 48CM/5MM/0.8S (**4/8 Slice**) from the test pull-down menu, and continue to slide cursor to the right to select AUTO.
 - d.) Click on the ACCEPT button.
 - 3.) Review the Image Analysis Report data and record the Series Means (AvXo-AvXc, and AvSDc) data and the Center Smudge data for each row in [Table 6-3](#) (**16 Slice**) or [Table 6-4](#) (**4/8 Slice**).
 - a.) Verify the Overall Test Pass/Fail Indicator Window of the Report display indicates PASS. (See [Table6-2](#).)
 - b.) Verify the scan parameter values displayed in the report match those in [Table6-2](#).
 - c.) Record the Series Means (AvXo-AvXc, and AvSDc) data and the Center Smudge data for each row in [Table 6-3](#) (**16 Slice**), or [Table6-4](#) (**4/8 Slice**).

6.4.2.2 Band and Streak Artifact

- 1.) Review the 48cm Phantom images acquired in Section 6.4.1 for Bands and Streaks using the Image Works Desktop.
 - a.) From the Global Control Palette, click on the IMAGE WORKS Desktop.
 - b.) From the Image Works Browser window, select the 48cm Phantom series acquired in Section 6.4.1.
 - c.) From the Image Works Browser window, click on the VIEWER button.
 - d.) While paging through the 64 images in the series, analyze images for any Band or Streak artifacts. Adjust Window/Level using the center mouse button, as required.
 - e.) Note any image that appears to have a Band or Streak artifact.
- 2.) Verify any images noted above (step 1e) meet 48cm Phantom Band or Streak Image Performance requirements.
 - a.) From the Global Control Palette, click on the Service Desktop.
 - b.) From the Service Browser window, select an image of the 48cm Phantom series acquired in Section 6.4.1 that appeared to have a Band or Streak while reviewing the series in Image Works.
 - c.) On the Image Analysis Tool window, click on the AUTO 1X rectangular button, select the IMAGE SERIES 48CM test from the pull-down menu, slide cursor to the right and select 16X1.25/120KV/200MA/2S (16 Slice) or 48CM/5MM/0.8S (4/8 Slice) from the test pull-down menu, and continue to slide cursor to the right to select BAND or STREAK.
 - d.) Click on the ACCEPT button.

Note: If required, adjust Window/Level using the center mouse button.

- e.) Position and size the Band or Streak ROI using the left cursor button.

Note: To ensure specifications are properly applied, do not adjust the Reference ROI's for the Band and Streak artifact tests. The Reference ROI's are adjustable after the 1st Accept Modification click.

- f.) Click on ACCEPT MODIFICATION twice, to generate a report.
- g.) Verify the Overall Test Pass/Fail Indicator Window of the Report display indicates PASS.
- h.) Repeat steps b through g, for each image noted in Step 1e.

Box Size =1600 mm2 to 2500 mm2

45 mm (+/- 5 mm) x 45 mm (+/- 5 mm)

45 (+/- 4 pixels) x 45 (+/- 4 pixels)

Box Positions:Box 1= 0 mm x 0 mm

Box 2 = 0 mm x -190 mm

Box 3 = 190 mm x 0 mm

Box 4 = 0 mm x 190 mm

Box 5 = -190 mm x 0 mm

Row	Images	Brightness Uniformity (AvXo - AvXc)	Noise (AvSDo)	Center Smudge Row Values	Comments
8A	1, 17, 33, 49				
7A	2, 18, 34, 50				
6A	3, 19, 35, 51				
5A	4, 20, 36, 52				
Specifications		< +/- 17	< 70	< +/- 14	

Table 6-3 48cm Phantom CT# Brightness Uniformity, Noise & Center Smudge Row Performance Data Sheet (16 Slice)

Row	Images	Brightness Uniformity (AvXo - AvXc)	Noise (AvSDo)	Center Smudge Row Values	Comments
4A	5, 21, 37, 53				
3A	6, 22, 38, 54				
2A	7, 23, 39, 55				
1A	8, 24, 40, 56				
1B	9, 25, 41, 57				
2B	10, 26, 42, 58				
3B	11, 27, 43, 59				
4B	12, 28, 44, 60				
5B	13, 29, 45, 61				
6B	14, 30, 46, 62				
7B	15, 31, 47, 63				
8B	16, 32, 48, 64				
Specifications		< +/- 17	< 70	< +/- 14	

Table 6-3 48cm Phantom CT# Brightness Uniformity, Noise & Center Smudge Row Performance Data Sheet (16 Slice)

Row	Images	Brightness Uniformity (AvXo - AvXc)	Noise (AvSDo)	Center Smudge Row Values	Comments
2A	1, 5, 9, 13				
1A	2, 6, 10, 14				
1B	3, 7, 11, 15				
2B	4, 8, 12, 16				
Specifications		< +/- 17	< 70	< +/- 14	

Table 6-4 48cm Phantom CT# Brightness Uniformity, Noise & Center Smudge Row Performance Data Sheet (4/8 Slice)

6.4.3 Failure Recommended Actions

6.4.3.1 Series Means (AvXo - AvXc or AvSDo) Failure Recovery

Specifications

Each Row of the series must pass 48cm Series Means (for the first series scan parameters) specifications:

- AvXo - AvXc: < +/- 17
- AvSDo: < 70.0

Recommended Recovery

- 1.) Repeat Detailed Cal for the 48cm Phantom.
- 2.) Repeat Sections 6.4.1 & 6.4.2, to verify Image Performance.

6.4.3.2 Center Smudge Failure Recovery

Specifications

80% the Rows must pass Center Smudge specifications when re-scanning:
Center Smudge:Smudge Factor < +/- 14.0

Recommended Recovery

- 1.) If one image of one Row fails (Row fails once):
 - Re-scan the 48cm Phantom acquiring two series of 64 images using the same protocol settings as listed in [Table 6-2](#).
 - Analyze the two 64 re-scanned image series using the Image Analysis IMAGE SERIES 48CM -> 16X1.25/120KV/200MA (**16 Slice**) or 48CM/5MM/0.8S (**4/8 Slice**) -> AUTO tool.
 - Verify that all the 128 re-scanned images do not fail Series Means (AvXo - AvXc and AvSDo) or Center Smudge.
- 2.) If more than one image of one Row fails (Row fails more than once) or multiple rows are failing:
 - Repeat Detailed Cal for the 48cm Phantom.
 - Re-scan the 48cm Phantom acquiring three series of 64 images using the same protocol settings as listed in [Table 6-2](#).
 - Analyze the three re-scanned image series using the Image Analysis IMAGE SERIES 48CM -> 16X1.25/120KV/200MA (**16 Slice**) or 48CM/5MM/0.8S (**4/8 Slice**) -> AUTO tool.
 - Verify that all the 192 re-scanned images do not fail Series Means (AvXo - AvXc and AvSDo) or Center Smudge.

6.4.3.3 Band or Streak Artifact Failure Recovery

Specifications

- Band Artifact: Band Factor< +/- 8.0
- Streak Artifact: Streak Factor< +/- 4.0

Recommended Recovery

- 1.) Troubleshoot Band or Streak artifact failures, as described in the *System Service Manual*.
- 2.) Repeat Detailed Cal for the 48cm Phantom.
- 3.) Repeat Sections 6.4.1 & 6.4.2, to verify Image Performance.

6.5 20cm QA Phantom Image Series Image Performance Verification

6.5.1 Acquiring the 20cm QA Phantom Image Series

- 1.) Mount the Phantom Holder on the head-end of the table.
- 2.) Mount the 20cm QA Phantom on the Phantom Holder.
- 3.) Align, level, & center the 20cm QA Phantom.
 - Align black line on phantom using the internal laser lights.
 - Level phantom using bubble level and the Z Axis knob on the Phantom Holder.

- Center phantom using the CENTER PHANTOM procedure in the left head SCANNER UTILITIES selection and the X and Y Axis knobs on the Phantom Holder.
- 4.) Set up the system to scan three 20cm QA Phantom image series with a Recon of the 1st series.
 - a.) On the Exam Rx desktop, select NEW PATIENT.
 - b.) Type the following entries in the listed Patient Information following fields:
 - * Patient ID: **Service**
 - * Name: **20cm QA Phantom Image Series**
 - c.) From the Protocol Selection display, click on the Service Tab under Anatomical Selector.
 - d.) On the Service Protocol Selection window, click on MANUFACTURING.
 - e.) On the Service Category Protocol List, click on the appropriate Image Series selection as follows:

System	Protocol - First	Protocol - Second
BrightSpeed Select	Image Series QA	(none)

- f.) On the ExamRx protocol parameter display, select the 1st Series.
- g.) Set internal Landmark.
- 5.) Acquire the 1st 20cm QA Phantom image series and one Recon series (MTF and Contrast, Visible Lines).
- 6.) Acquire the 3rd 20cm QA Phantom image series.
- 7.) Acquire the last 20cm QA Phantom image series.

6.5.2 20cm QA Phantom Image Series Image Performance Verification

1st QA Phantom Image Series (4 Image MTF Average and Contrast Scale) Image Performance Verification

- 1.) Launch the Image Analysis Tool.
 - a.) From the Global Control Palette, click on the Service Desktop.
 - b.) From the Service Desktop upper navigational bar, click on IMAGE QUALITY.
 - c.) In the left navigational area directory, click on IMAGE ANALYSIS.
- Note: During the Image Analysis process, you will need to switch between the Service Browser screen (to select the image series to analyze) and the Image Analysis screen (to run a particular Auto 1x analysis test). Use ALT/F1 to bring a window to the front. Use ALT/F3 to send a window to the back.
- 2.) Using the Auto 1x Tool on the Image Analysis screen, generate a QA#1 MTF and QA#1 Contrast Scale report for the 1st series of the 20cm QA Phantom acquisition.
 - a.) From the Service Browser window, select the 1st series of the 20cm QA Phantom series acquired in section 6.5.1.
 - b.) On the Image Analysis Tool window, click on the AUTO 1X diamond shaped button. (See [Figure 6-2, on page 234.](#))
 - c.) Click on the AUTO 1X rectangular button, select the IMGSER 20QA test from the pull-down menu, slide cursor to the right to select QA#1 from the test pull-down menu.
 - d.) Click on the ACCEPT button.
 - 3.) Review the Image Analysis Report data and record the Per Image MTF, 4 Image MTF Average, and Per Image Contrast Scale data for each row in [Table 6-5](#).
 - a.) Verify the Overall Test Pass/Fail Indicator Window of the Report display indicates PASS.

(See [Figure 6-3, on page 235.](#))

- b.) Record Per Image MTF and Per Image Contrast Scale data for each image in [Table 6-5](#).
- c.) Record the 4 Image MTF Average in [Table 6-5](#).

Image	MTF	MTF 4-slice average	Contrast Scale	Comments
1		N.A.		
2		N.A.		
3		N.A.		
4		N.A.		
Specifications	N.A.	0.65 to 1.0	110.0 to 130.0	N.A.

Table 6-5 20cm QA#1 Phantom High Contrast Spatial Resolution Image Performance (MTF and Contrast Scale)

2nd QA Phantom Image Series (4 Bone Retro Image Visible Lines) Image Perf. Verification

- 1.) Using the Auto 1x Tool on the Image Analysis screen, generate a Per Image QA#1 Lines report for the 2nd series of the 20cm QA Phantom acquisition.
 - a.) From the Service Browser window, select the 1st image of the 2nd series of the 20cm QA Phantom series acquired in section 6.5.1 (on [page 239](#)).
 - b.) Click on the AUTO 1X rectangular button, select the IMGSER 20QA test from the pull-down menu, and slide cursor to the right to select QA#1/LINES/0.8S.
 - c.) Click on the ACCEPT button. The tool displays a Visible Lines pop-up window.
 - d.) On the Visible Lines pop-up window, click on the VISIBLE LINE button.
- Note: The Image Analysis Tool automatically sets the Window/Level for optimal viewing. If required, adjust Window/Level using the center mouse button.
- e.) From the Visible Line pull-down menu, select the letter (A through F) that matches the smallest line pair pattern that you can discern in the Image View Port.
- f.) On the Visible Lines pop-up window, click on the OK button.
- 2.) Review the Image Analysis Report data and record the Per Image QA#1 Visible Lines data for the 1st image in [Table 6-6](#).
 - a.) Verify the Overall Test Pass/Fail Indicator Window of the Report display indicates PASS. (Refer to [Figure 6-3, on page 235.](#))
 - b.) Record the Per Image QA#1 Visible Lines data for the 1st image in [Table 6-6](#).
- 3.) Repeat Steps 1 and 2, for each of the remaining three images in the 2nd 20cm QA Image Series.

Image	Line Patterns Visible	Comments
1		
2		
3		
4		
Specifications	B, C, D, E, F	N.A.

Table 6-6 20cm QA#1 Phantom High Contrast Spatial Resolution Image Performance (Visible Lines)

3rd QA Phantom Image Series (4 Image Visible Holes) Image Performance Verification

- 1.) Using the Auto 1x Tool on the Image Analysis screen, generate a Per Image QA#2 Holes report for the 3rd series of the 20cm QA Phantom acquisition.
 - a.) From the Service Browser window, select the 1st image of the 3rd series of the 20cm QA Phantom series acquired in section 6.5.1.
 - b.) Click on the AUTO 1X rectangular button, select IMGSER 20QA from the pull-down menu, and slide cursor to the right to select QA#2 HOLES.
 - c.) Click on the ACCEPT button. The tool displays a Visible Hole pop-up window.
 - d.) On the Visible Hole pop-up window, click on the VISIBLE HOLE button.
- Note: The Image Analysis Tool automatically sets the Window/Level for optimal viewing. If required, adjust Window/Level using the center mouse button.
- e.) From the Visible Hole pull-down menu, select the number (1 through 5) that matches the smallest hole that you can discern in the pattern centered in the image displayed in the Image View Port.
- f.) On the Visible Hole pop-up window, click on the OK button.
- 2.) Review the Image Analysis Report data and record the Per Image QA#2 Visible Holes data for the 1st image in [Table 6-7](#).
 - a.) Verify the Overall Test Pass/Fail Indicator Window of the Report display indicates PASS. (Refer to [Figure 6-3](#), on [page 235](#).)
 - b.) Record the Per Image QA#2 Visible Holes data for the 1st image in [Table 6-7](#).
- 3.) Repeat Steps 1 and 2 for the 3rd, 5th, and 7th images in the 3rd 20cm QA Image Series.

Image	Visible Holes Viewable at Window 20	Contrast Factor	Comments
1			
3			
5			
7			
Specifications	See Table 6-9	2.0 to 12.0	N.A.

Table 6-7 20cm QA#2 Phantom Low Contrast Detectability Image Performance (Visible Holes)

Contrast Factor Range (Box 1 Means - Box 2 Means)	Visible Number of Holes		Smallest Visible Hole Size
	Lower Limit *	Upper Limit *	
2.00 to 3.99	2	5	7.5 mm
4.00 to 7.99	3	5	5.0 mm
8.00 to 12.00	4	5	3.0 mm

* Required Number of Visible Holes depends on the Contrast Factor

Table 6-8 20cm QA#2 Phantom Visible Hole Specifications

Last QA Phantom Image Series (QA#3 Small) Image Performance Verification

- 1.) Using the Auto 1x Tool on the Image Analysis screen, generate a QA#3 Small report for the last series of the 20cm QA Phantom acquisition.
 - a.) From the Service Browser window, select the last series of the 20cm QA Phantom series acquired in section 6.5.1 (on [page 239](#)).

- b.) Click on the AUTO 1X rectangular button, select IMGSER 20QA from the pull-down menu, slide cursor to the right to select QA#3/SMALL/1.0S from the test pull-down menu, and continue to slide cursor to the right to select AUTO. (Refer to [Figure 6-3, on page 235](#).)
 - c.) Click on the ACCEPT button.
 - 2.) Review the Image Analysis Report data and record the QA3 Small data for the last 20cm QA Image Series in [Table 6-9](#).
 - a.) Verify the Overall Test Pass/Fail Indicator Window of the Report display indicates PASS. (Refer to [Figure 6-3, on page 235](#).)
 - b.) Record the QA3 Small data (Row 2A1A and Row 2B1B Series Means AvXc and AvXo - AvXc), QA3 (AvSDc), and the Center Smudge Row data in [Table 6-9](#).

Box Size =196 mm2 to 256 mm2

15 mm (+/- 1 mm) x 15 mm (+/- 1 mm)

31 (+/- 2 pixels) x 31 (+/- 2 pixels)

Box Positions:Box 1= 0 mm x 0 mm
 Box 2 = 0 mm x -80 mm
 Box 3 = 80 mm x 0 mm
 Box 4 = 0 mm x 80 mm
 Box 5 = -80 mm x 0 mm

Row	Images	AvXc	AvXo	AvXo - AvXc	AvSDo	AvSDc	Center Smudge Row Value	Comments
2A1A	1, 3, 5, 7				N.A.			
1B2B	2, 4, 6, 8				N.A.			
Specifications		+3.0 to -3.0	N.A.	+3.0 to -3.0	N.A.	+ 2.9 to +3.5	+2.2 to -2.2	

Table 6-9 20cm QA#3 Phantom CT#, Brightness Uniformity, Center Noise & Center Smudge Row Performance Data Sheet

- 3.) Review the last 20cm QA Phantom images acquired in Section 6.5.1 for Rings, Bands, and Streaks using the Image Works Desktop.
 - a.) From the Global Control Palette, click on the IMAGE WORKS DESKTOP.
 - b.) From the Image Works Browser window, select the last 20cm QA Phantom series acquired in section 6.5.1.
 - c.) From the Image Works Browser window, click on the VIEWER button.
 - d.) While paging through the eight images in the series, analyze images for any Ring, Band, or Streak artifacts. Adjust Window/Level using the center mouse button, as required.
 - e.) Note any image that appears to have a Ring, Band, or Streak artifact.
 - 4.) Verify any images noted above (step 3e) meet 20cm QA Phantom Ring, Band, or Streak Image Performance requirements.
 - a.) From the Global Control Palette, click on the Service Desktop.
 - b.) From the Service Browser window, select an image of the 4th 20cm QA Phantom series acquired in section 6.5.1 that appeared to have a Ring, Band, or Streak while reviewing the series in Image Works.
 - c.) On the Image Analysis Tool window, click on the AUTO 1X rectangular button, select the IMGSER 20QA test from the pull-down menu, slide cursor to the right to select QA#3/SMALL/1.0S from the test pull-down menu, and continue to slide cursor to the right to select RING, BAND, or STREAK.
 - d.) Click on the ACCEPT button.
- Note: If required, adjust Window/Level using the center mouse button.
- e.) Position and size the Ring, Band, or Streak ROI, using the left cursor button.

- Note: To ensure specifications are properly applied, do not adjust the Reference ROIs for the Ring, Band, and Streak artifact tests. The Reference ROIs are adjustable after the 1st Accept Modification click.
- f.) Click on ACCEPT MODIFICATION twice to generate a report.
 - g.) Verify the Overall Test Pass/Fail Indicator Window of the Report display indicates PASS.
 - h.) Repeat steps b through g for each image noted in Step 3e.

6.5.3 Failure Recommended Actions

6.5.3.1 1st Image Series (4 Image MTF Average) Failure Recovery

Specifications

The 4 Image MTF Average must pass specifications:

4 Image MTF Average: 0.65 to 1.00

Recommended Recovery

- 1.) Repeat Detailed Cal for the 20cm QA Phantom.
- 2.) Repeat Auto CT# Adjust for the 20cm QA Phantom.
- 3.) Repeat Sections 6.5.1 & 6.5.2, to verify Image Performance.

6.5.3.2 2nd Image Series (Visible Lines) Failure Recovery

Specifications

The Largest 5 line pair patterns (coded F, E, D, C, and B) must be visible for each of the four images in this series.

Recommended Recovery

Note: The most common failure for this test is that the phantom has air bubbles that are obscuring the line pair patterns.

- 1.) Carefully inspect the 20cm QA Phantom for air bubbles. If required, refill the phantom with water, to eliminate all air bubbles.
- 2.) Repeat Sections 6.5.1 & 6.5.2, to verify Image Performance.

6.5.3.3 3rd Image Series (Visible Holes) Failure Recovery

Specifications

The number of visible holes in this image series must pass specifications as follows:

- Contrast Factor = 2.0 to 3.99:
Minimum of 2 holes must be visible (10mm & 7.5mm)
- Contrast Factor = 4.0 to 7.99:
Minimum of 3 holes must be visible (10mm, 7.5mm & 5mm)
- Contrast Factor = 8.0 to 12.0:

Minimum of 4 holes must be visible (10mm, 7.5mm, 5mm, & 3mm)

Recommended Recovery

- 1.) Repeat Detailed Cal for the 20cm QA Phantom.
- 2.) Repeat Auto CT# Adjust for the 20cm QA Phantom.
- 3.) Repeat Sections 6.5.1 & 6.5.2, to verify Image Performance.

6.5.3.4 4th Image Series Failure Recovery

Series Means

SPECIFICATIONS

Both Rows (2A1A, 1B2B) of the series must pass QA#3 Series Means specifications:

- AvXc: < +/- 3.0
- AvXo - AvXc: < +/- 3.0

RECOMMENDED RECOVERY

- 1.) Repeat Detailed Cal for the 20cm QA Phantom.
- 2.) Repeat Auto CT# Adjust for the 20cm QA Phantom.
- 3.) Repeat Sections 6.5.1 & 6.5.2, to verify Image Performance.

QA#3 (AvSDc)

SPECIFICATIONS

Both Rows (2A1A, 1B2B) of the series must pass QA#3 AvSDc specifications:

- AvSDc (< 5K Scans):2.90 to 3.50
- AvSDc (> 5K Scans):2.80 to 3.60

RECOMMENDED RECOVERY

- 1.) Repeat Detailed Cal for the 20cm QA Phantom.
- 2.) Repeat Auto CT# Adjust for the 20cm QA Phantom.
- 3.) Repeat Sections 6.5.1 & 6.5.2, to verify Image Performance.

Center Artifact or Center Smudge

SPECIFICATIONS

All eight images of the series must pass 20cm QA#3 Center Artifact and Center Smudge specifications:

- Center Artifact Factor: +/- 3.5
- Center Smudge Factor: +/- 2.2

RECOMMENDED RECOVERY

- 1.) Repeat Detailed Cal for the 20cm QA Phantom.
- 2.) Repeat Auto CT# Adjust for the 20cm QA Phantom.
- 3.) Repeat Sections 6.5.1 & 6.5.2, to verify Image Performance.

Ring, Band, or Streak Artifact Failure Recovery

SPECIFICATIONS

- Ring Artifact: Ring Factor < +/- 4.8
- Band Artifact: Band Factor +/- 2.8 from 0cm to 8.5cm
< +/- 2.8 > 8.5cm
- Streak Artifact: Streak Factor < +/- 4.0

RECOMMENDED RECOVERY

- 1.) Troubleshoot Ring, Band or Streak artifact failures as described in the Service Manual.
- 2.) Repeat Detailed Cal for the 20cm QA Phantom.
- 3.) Repeat Sections 6.5.1 & 6.5.2, to verify Image Performance.

Section 7.0

System Functional Test

Use the system tests in the following sections to exercise all aspects of the system and to ensure system integrity before releasing to the customer. Although the means, standard deviation, and resolution specifications do not apply during system functional tests, treat any artifact or image anomaly as a failure.

If you encounter a failure during the system tests:

- Record any evidence of artifacts, such as rings, streaks, shading, cupping, noise, or center artifacts.
 - Correct artifacts, system test, or image series failures when they occur. Any delay in repairs could increase the number of retests.
- 1.) Place the QA phantom on the cradle.
 - Drive the table to an elevation of 100.
 - Align the line on the phantom with the internal laser lights.

Never scan above 50mA without first placing a phantom in the field of view. Levels in excess of 50mA can cause temporary radiation damage to the detector that lasts several hours. If you acquire image series cals with a radiation damaged detector, the cals may cause artifacts in subsequent image series scans.

- 2.) Select the service protocol, SYSTEM SCAN/CUST QUAL RELI.
- 3.) Run each series.

Note:

Stop the service protocol after the second helical series. It is not necessary to proceed past that point in the protocol.

- 4.) Using Diagnostic Data Analysis, review the MSD Plots for each AIR scan.
- 5.) Complete the scans.

Section 8.0

Save System State

Use the following commands to create the System State DVD.

- 1.) Load a DVD into the mod drive on the front of the console.

- 2.) If you are not on the Service Desktop, click on the SERVICE DESKTOP icon, .

- 3.) Click on the UTILITIES icon, .

- 4.) Select SYSTEM STATE to open the System State Save/Restore menu.

- 5.) Select ALL

- 6.) Select SAVE

- 7.) When the save operation completes, select FILE and QUIT from the pull down menu.

- 8.) Remove the DVD from the drive

Chapter 7

Customer Options Installation & Verification

Note: Only use the Installation manual that arrives with your system for installation. Any other revisions of this manual may not exactly match your system.

Section 1.0

CT Options

- ☐ Injector functional tests completed.
- ☐ AWW functional tests completed.
- ☐ Filming/Camera/DASM functional tests completed.
- ☐ UPS functional tests completed.
- ☐ Network items installed and functional tests completed.
- ☐ Verify that all customer software options are installed and functional.

Section 2.0

Install Options

Refer to:

- GE Prints and schematics for mechanical (physical) location of option
- FDO shipment for identification of items
- Installation Specialist for installation instructions if they differ from print

2.1 5 inch MOD

Refer to documentation shipped with 5" MOD.

- 5" MOD Drive Option 5162230

2.2 Camera (Filming Device)

Refer to documentation shipped with camera

- Chapter 2 of this manual
- DICOM 2210573 GE Document
- DICOM Print 2152913

2.3 Advantage Windows Workstation (AWW)

Refer to the directions provided:

- Pre-install 2111833
- Service 2111831

Section 3.0

DICOM Network Introduction

BrightSpeed Select systems support two basic Networking Protocols:

- Advantage NET (IC, Signa 4.X, CT-HLA, CT/I ...)
- DICOM (CT/I, CT Synergy, Advantage Workstations, ...)

DICOM NETWORK INTRODUCTION

DICOM networks basically operate on the **tasks** or services that various devices on the network use or provide. These services are labeled as Application Entity Titles (AE Titles). The CT scanner system is a **user** of six DICOM Network Services and is a **provider** of two DICOM Services:

AS A DICOM SERVICE USER:

- Send or **Push** images to another network device.
- Send or **Push** images to a DICOM Printer.
- Review image database on another device and retrieve or **Pull** selected images from that device (Query/Retrieve User).
- Send or **Push** images to a an image storage device and obtain confirmation that the images have been archived (Storage Commitment).
- Obtain Patient Worklist Information from the Hospital HIS/RIS System.
- Store images on DVD-RAM media.

AS A DICOM SERVICE PROVIDER:

- Receive **Pushed** images from another network device
- Allow another network device to review the image database and to retrieve or **Pull** selected images (Query/Retrieve Provider)

For each DICOM Service that the CT system will be a **User** (except for storing images on DVD-RAM media), you must **declare** this device on the CT system using three menu selections. For some devices, you must declare not only the device, but each service (AE Title) that the device provides. For example, you may be required to declare a PACS System twice on the CT system: once as a destination to **push** images and, second, as destination that provides storage commitment capability after images have been **pushed**.

For each DICOM Service that the CT system will be a **Provider**, you must **declare** the CT system on the network device that will be using these services.

Information required to complete configuring a hospital DICOM network is provided by the hospital network administrator (hostnames, IP Addresses) and the DICOM Conformance Statement document (AE Titles, Port Numbers) provided with each DICOM compatible network device on the network.

Section 4.0

Before You Start

4.1 Network Physical Requirements

Before setting up the CT scanner system on the hospital network, verify the following physical items are complete:

- Scanner console, monitor, keyboard, and mouse are installed and connected.
- CT system power is ON
- Hospital Ethernet network RJ45 Class IV twisted pair cable is connected to the scanner console network receptacle.
- Hospital network connection is operational and is running 10baseT or 100baseT.

4.2 Network Identity Information

To declare the CT system on the network, ensure the following network identity information is available from the Hospital Network Administrator:

- Hostname (No more than 16 Characters)
- Internet Protocol (IP) Address
- Subnet Net Mask IP Address (if applicable)
- Broadcast Address (if applicable)

4.3 Scanner to DICOM Remote Hosts Network Information

To declare DICOM remote hosts (PACS systems, archival devices, review workstations) on the CT system, ensure the following information is available for each remote host:

From the Hospital Network Administrator:

- Hostname
- Internet Protocol (IP) Address
- Network Protocol (DICOM for CT Systems)

From the Remote Host Device DICOM Conformance Statement Document:

- DICOM Application Entity Title or AE Title (DICOM service that remote host provides or uses)
- DICOM Listening Port Number

4.4 Scanner to DICOM HIS/RIS Interface Network Information

To declare DICOM Hospital HIS/RIS Interface devices (Mitra and others) on the CT System, ensure the following information is available:

From the Hospital Network Administrator: *Internet Protocol (IP) Address*

From the HIS/RIS Interface Device DICOM Conformance Statement Document:

- DICOM Application Entity Title or AE Title (DICOM Service that the HIS/RIS interface provides)
- DICOM Listening Port Number

4.5 Scanner to DICOM Printers Network Information

To declare DICOM on the CT System, ensure the following information is available for each printer:

From the Hospital Network Administrator:

- Hostname
- Internet Protocol (IP) Address

From the Printer DICOM Conformance Statement Document:

- DICOM Application Entity Title or AE Title (DICOM service that remote host provides or uses)
- DICOM Listening Port Number

Section 5.0

Declaring the System on the Hospital Network

5.1 Enter Configuration Routine

- 1.) On the OC, open a Unix shell window.
 - 2.) Enter root as a superuser:
Type: **su -** ENTER at the prompt.
Type: **#bigguy** ENTER at the password prompt.
 - 3.) Change directory to scripts:
Type: **cd /user/g/scripts** ENTER at the root prompt.
 - 4.) Launch the Install Utility:
Type: **reconfig** ENTER at the prompt.
- The OC displays the Install Utility Window as shown in [Figure 7-1](#).

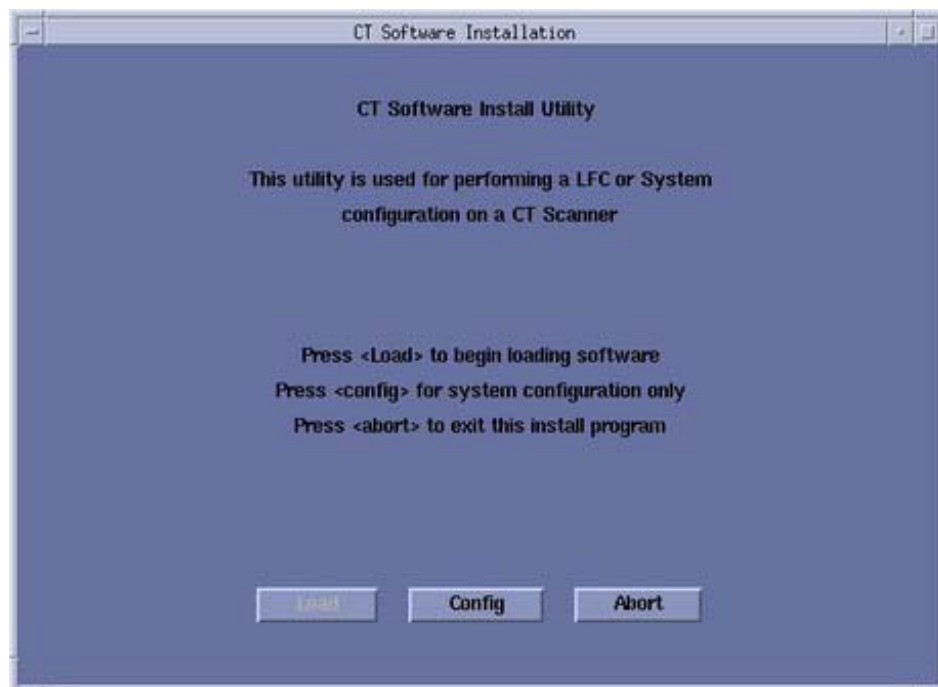


Figure 7-1 Install Utility Window

5.) Enter the Configuration Routine:

Using the mouse, click on the CONFIG button.

The OC displays the System Configuration - System Settings screen, as shown in [Figure 7-2](#).

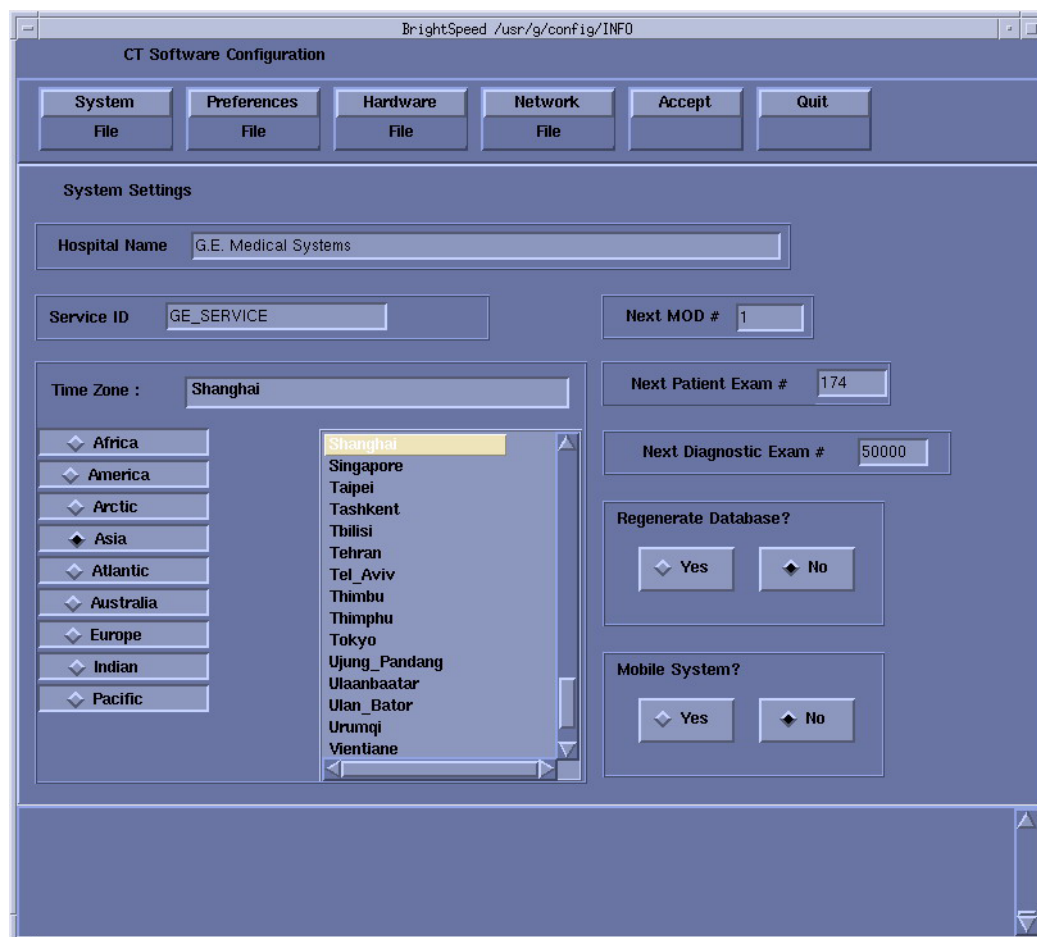


Figure 7-2 System Settings Screen

5.2 Configure Network Settings

This screen provides the ability to declare the CT system on a hospital network. Key information such as Host Name, IP Address, Net Mask (for CT systems on a subnet) must be obtained from the hospital network administrator.

- 1.) Select the NETWORK button to display the Network Settings screen as shown in [Figure 7-3](#).
- 2.) Enter the Suite Name.

The Suite Name is a means of identifying this particular CT system as a part of a group of CT Systems in a suite configuration. This Suite Name will appear on all image headers.

The Suite Name must start with a letter, followed by three alphanumeric characters (total MUST be four characters long). The name of the OC interface will be <Suite Name>_oc and the SBC interface will be <Suite Name>_sbc.

CT Software Configuration

System Preferences Hardware Network Accept Quit

File File File File

Network Settings

Gateway Parameters

Suite Name CT99

Host Name ct99

IP Address 3.36.231.128

Net Mask 255.255.252.0

Broadcast Address 3.36.231.255

Default Gateway 3.36.231.1

☐ Advanced Options

Figure 7-3 Networks Settings Screen

- 3.) Enter the hospital provided Host Name.
The Host Name identifies the network hostname and AE Title of the CT system.
The Host Name:
 - **MUST NOT** be <Suite Name>_oc or <SUITE NAME>_OC.
 - **MUST NOT** exceed 16 Characters.
 - **MUST** only contain the following characters: **A** through **Z**, **a** through **z**, **0** through **9**, - and _
- 4.) Enter the hospital provided IP Address.
- 5.) Enter the hospital provided Net Mask (if the CT system is on a subnet).
- 6.) Enter the Broadcast Address:
The Broadcast Address should be the same as the IP Address except for the bits of the host id portion (last digit group) set to 1s or 0s depending on the configuration of the network. The standard default is 1s but older SunOS machines used 0s.
Example: If the IP Address is 192.100.9.17, the Broadcast Address should be 192.100.9.255 if the network is configured to use 1's to specify the broadcast address.
If the network contains genesis based scanners or other SunOS 3.5 or 4.1 computers, the Broadcast Address should be 192.100.9.0.
- 7.) Enter the hospital provided Default Gateway IP Address in the Default Gateway field (if applicable). If the site network does not use a default gateway, leave the field blank.

- 8.) Select NIS (Yellow Pages database) Advanced Option only if requested by the hospital network administrator as follows:
 - a.) Select ADVANCED OPTIONS button on the Network Settings screen.
 - b.) Select Use NIS? button.
 - c.) Enter the hospital provided Domain Name.
- 9.) Record all the Network parameters in the *Software Installation Procedures* Document, or on the worksheet in 2.2 on page 199.

5.3 Initiate System Reconfiguration

- 1.) Select ACCEPT on the System Configuration Screen.

The system loads the application software, OS patches, and kernel changes, and configures the system on both the OC and the SBC.

This loading process takes approximately 15 minutes. While the load is going on, the results are displayed in a shell window, which closes when the loading process is complete. All the window output is logged to a file named:

`/var/adm/install.log.YYYYMMDDWWHHMMSS.`

(Where `YYYYMMDDWWHHMMSS` is the Date/Time that the loading process was started.)
- 2.) When the loading process and configuration changes are complete, the system displays a prompt to reboot. Click on YES.
- 3.) The system will automatically login as ctuser after the reboot. Select OK on the Autostart Disabled popup message.
- 4.) To startup Applications, in the console shell window, type **startup** ENTER.

Section 6.0

Declaring Remote Hosts on the CT System

6.1 Enter Remote Host Configuration Screen

- 1.) On the OC, select the IMAGE WORKS icon.
- 2.) Select NETWORK.



6.2 Declaring Advantage NET Remote Hosts on the Scanner

Use Advantage NET Protocol networks to communicate with older CT or MR Systems (MR Signa 4.x, CT-HLA, CT/I Systems, and Workstations that support the Advantage NET protocol). Advantage NET Protocol does not offer full compatibility with BrightSpeed Select DICOM formats.

Repeat the following procedure for each Advantage NET Remote Host device that the customer expects to have this CT system communicating with.

- 1.) Select REMOTE HOSTS from the pull down menu. The system displays the Remote Host Parameter Screen as shown in [Figure 7-4](#).

Figure 7-4 Advantage Net Network Protocol Parameter Settings

- 2.) Enter the hospital provided Host name.
- 3.) Enter the hospital provided Network Address (IP Address).

- 4.) Select ADVANTAGE NET as the Network Protocol.

The systems automatically de-highlights the remaining parameter fields on the Remote Host parameter selection screen. These are dedicated DICOM protocol parameters and do not apply to Advantage NET type devices.

- 5.) Select SAVE to store the parameter settings of the remote host.

6.3 Declaring DICOM Remote Hosts on the CT Scanner

Use DICOM protocol networks to communicate to DICOM devices such as CT/i, CT Synergy, DLX, MR Lx, and third party hosts.

Repeat the following procedure for each DICOM remote host device that the customer expects to have this CT system communicating with.

- 1.) Select REMOTE HOSTS from the pull down menu. The system displays the Remote Host Parameter screen as shown in [Figure 7-5](#).

Figure 7-5 DICOM Network Setting Protocol Parameter Settings

- 2.) Enter the hospital provided Host name.
- 3.) Enter the hospital provided Network Address (IP Address).
- 4.) Select DICOM as the Network Protocol.
The system automatically highlights the remaining parameter fields on the Remote Host parameter selection screen. These are dedicated DICOM Protocol parameters.
- 5.) Enter the TCP/IP Listening Port Number from the DICOM Conformance Statement provided with the device.

- 6.) Enter the AE Title from the DICOM Conformance Statement provided with the device.
Application Entity Titles (also known as ACR-Nema or Dicom Name) refer to the DICOM Network Services that a device provides to the CT System. For most devices, the AE Title is the same as the hostname (CT systems are equipped with this feature).
However, some devices such as PACS systems may have separate AE Titles and port numbers for each of the services that the PACS system provides. In these cases, you must enter a separate remote host (same hostname and IP Address) for each of the independent AE Title Services that the host provides (one host as an image **push-to** destination, another host as a **query/retrieve** provider, and another host as a **storage/commitment** provider).
Be sure to review the DICOM Conformance Statement for each device that will provide a remote host network service for the CT system (image **push-to** or store destination, Query/Retrieve, and Storage Commitment) to ensure that each service is correctly configured.
- 7.) Select the correct Archive Node choice for the device. The Archive Node selection field defines the ability of the remote host to act as a DICOM Storage/Commitment provider and indicate to the operator that a study/series/image was archived. Select:
 - AUTO to have the CT system automatically check to see if the designated remote host is a DICOM Storage/Commitment Provider.
 - YES if the device is the hospital designated DICOM Storage/Commitment Provider. During an Application Study Archive process, the local browser screen will indicate Archive Status = Y to the operator.
 - NO if the device is not a DICOM Storage/Commitment Provider.
- 8.) Select the correct Access to local host: settings. These two selections allow you to selectively block the remote host from using the BrightSpeed Select DICOM services as a provider (image **push-to** destination, and a **Query/Retrieve** provider).
 - Send Images? Set to YES if the customer wants the CT system to be able to have images **pushed** to the system from the applicable remote host. Set to NO if the customer wants to block an image **push** from the applicable remote host.
 - Query/retrieve images? Set to YES if the customer wants the remote host to be able to review the image database (query) and **pull** selected images from the database. Set to NO if the customer does not want the remote host to have this ability.
- 9.) Select the correct Custom search? setting. This selection allows the CT scanner to selectively search through the remote host's image database when the operator is using remote browser screen to **query** the remote host. The search parameters that the CT system allows the customer to use are: last name contains, patient ID, exam number, accession number, and exam date.
 - Select ON if the device supports custom searches as part of the devices **Query/Retrieve** DICOM Provider service.
 - Select OFF if the device does not support custom searches.
- 10.) Record all the remote host network parameters for each remote host in the *Software Installation Procedures* Document.
- 11.) Select SAVE to store the parameter settings of the remote host.

Section 7.0

Declaring the CT System on Remote Hosts

7.1 Declaring the Scanner on Advantage NET Protocol Devices/Systems

Refer to the appropriate service manual provided with the Advantage NET Protocol device or system to find instructions how to declare the CT System as an Advantage NET remote host.

7.2 Declaring the Scanner on DICOM Protocol Devices/Systems

Refer to the appropriate Service Manual provided with the DICOM protocol device or system to find instructions how to declare the CT System as a DICOM remote host.

The CT System provides two DICOM Services as a provider to remote hosts:

- A remote host can **push** images to the CT image database.
- A remote host can review the CT image database (query) and **pull** selected images (retrieve).

Use the following parameter information to configure the DICOM device/system to either **push** images to the CT scanner and/or perform a **Query/Retrieve** operation:

- **Hostname**: Provided by the Hospital Network Administrator. Exactly the same scanner assigned hostname entered in Network Configuration Screen.
- **Application Entity Title**: Exactly the same entry as the Hostname.
- **Network Address**: Provided by the Hospital Network Administrator. Exactly the same scanner assigned IP Address entered in Network Configuration Screen.
- **Network Protocol**: DICOM 3.0.
- **Port Number**: For all DICOM service that the CT System provides, use 4006.
- **Provider Type**: This field concerns the BrightSpeed Select DICOM **Query/Retrieve** provider capability. All CT systems are wstudy root systems, which allow queries at the exam, series, and image level.
- **Support Worklist**: This field concerns whether a DICOM **Query/Retrieve** provider capable device or system supports a filter search of the image database. All CT systems support a filtered search of the image database as part of the BrightSpeed Select DICOM **Query/Retrieve** provider capability.

Section 8.0 DICOM HIS/RIS Setup

8.1 Prerequisites

Most hospital HIS/RIS systems are not DICOM compatible and require a DICOM HIS/RIS Worklist Interface to provide patient scheduling information to the CT system. Contact your local HNS support engineer to determine exactly what DICOM HIS/RIS Interface is appropriate for the customer.

In addition, the CT system must have the ConnectPRO software option installed to utilize the DICOM Protocol Worklist capability.

8.2 Loading ConnectPRO Software Option on the CT System

- 1.) If you are not on the Service Desktop, click on the SERVICE DESKTOP icon.



- 2.) Click on the CONFIGURATION icon.



- 3.) Click on INSTALL OPTIONS.

- 4.) Select INSTALL and click START. The console displays the Software Options window as shown in [Figure 7-6](#).

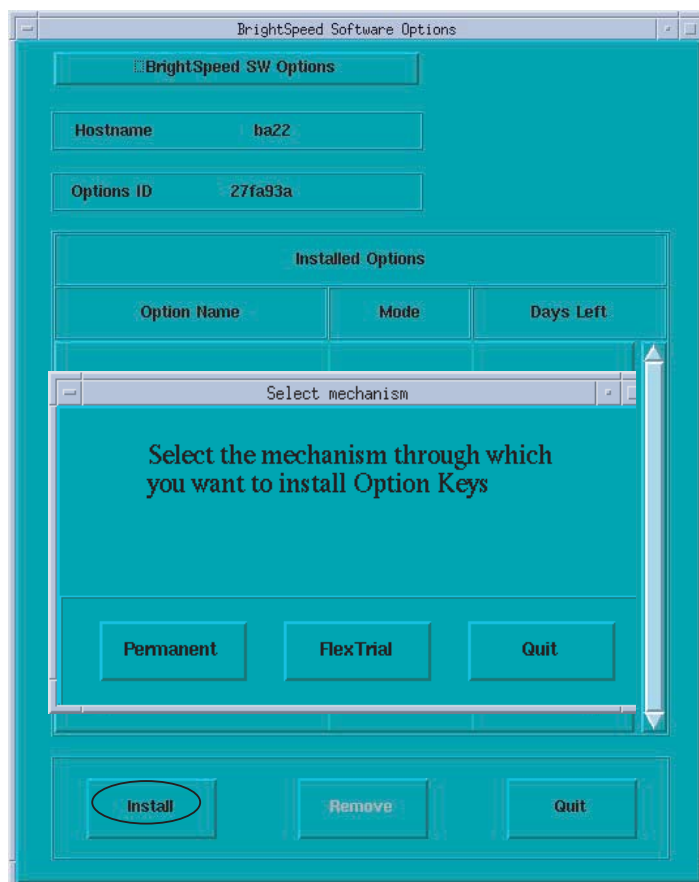


Figure 7-6 Options Window when First Selected

- 5.) Insert the ConnectPRO Options DVD into the DVD drive.
- 6.) Click on OK to continue.
- 7.) On the Software Option window, select the CONNECTPRO option in the Available Options window.
- 8.) Click on the INSTALL button.

The console may display a message box during the software loading operation. When the system has completed loading the software:

- the Installed Options window displays the ConnectPRO option, and
- the console displays the ConnectPro Setup window as shown in [Figure 7-7](#).

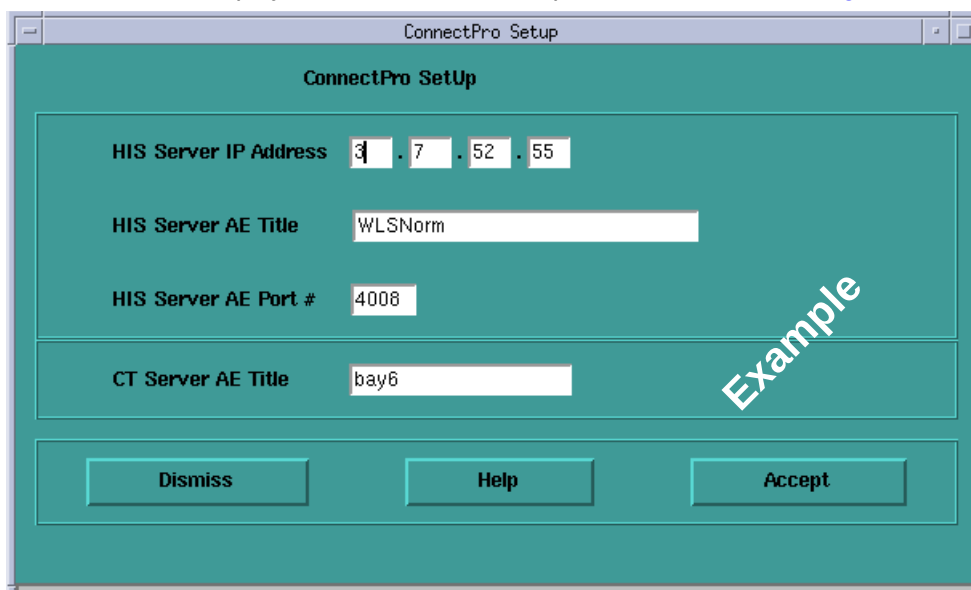


Figure 7-7 ConnectPro Setup Screen

- 9.) Enter the hospital provided HIS Server IP Address.
- 10.) Enter the HIS Server AE Title from the DICOM HIS/RIS Interface device DICOM Conformance Statement document.
- 11.) Enter the HIS Server AE Port # from the DICOM HIS/RIS Interface device DICOM Conformance Statement document.
- 12.) Record all the ConnectPRO HIS/RIS Network parameters.
The system will automatically load the BrightSpeed Select CT Server AE Title (Hostname of the system). CT Server AE Title will be used for PPS setup later.
- 13.) Click on the ACCEPT button.
- 14.) "Should PPS be enabled?" is popped up. Select Yes.
- 15.) The console displays the ConnectPro Setup window as shown in

PPS Server Setup

PPS Server (SCP) SetUp

PPS Server IP Address 3 . 70 . 204 . 86

PPS Server AE Title DICOM_TEST_RMT

PPS Server AE Port # 4500

CT Server AE Title bay6

Dismiss Help Accept

Example

Figure 7-8 PPS Setup Screen

- 16.) Enter the hospital provided PPS Server IP Address.
- 17.) Enter the PPS Server AE Title from the DICOM HIS/RIS Interface device DICOM Conformance Statement document.
- 18.) Enter the PPS Server AE Port # from the DICOM HIS/RIS Interface device DICOM Conformance Statement document.
- 19.) Enter CT Server AE Title got from [Step 12](#)). Click [Accept].
- 20.) The system will check valid setup entries, execute a perl script that makes changes to the configuration files for HIS/RIS, and return to the Software Option window shown in [Figure 7-6](#).
- 21.) Click on the QUIT button and the subsequent OK button in the message window.
- 22.) Remove the DVD from the DVD drive, and write protect the side containing the ConnectPRO option.

8.3 PPS Setup

If PPS is not enabled during ConnectPro Option installation, PPS can also be setup when needed later.

From the tool chest located in upper right corner of either monitor, select Unix Shell.

Type '**hostname**'. The output is the CT scanner AE title that will be sent to the worklist server.

Type '**installhisris**'.

Then the ConnectPro Setup window as shown in [Figure 7-7](#). Follow the [Step 9\)](#) to [Step 19\)](#) in [Section 8.2 on page 262](#).

Section 9.0

DICOM Filming Devices Setup

9.1 Prerequisites

Before configuring DICOM filming devices (cameras, printers) on the CT System, ensure the following are complete:

- Filming Device Service Representative to assist in camera/printer setup for best image quality film presentation.
- Hospital DICOM network is operational.
- Filming device is connected to the DICOM network with the correct filmer DICOM interface.
- Filming device is DICOM protocol compatible.
- Filming device DICOM Conformance Statement document is available.

Note: **Filmer DICOM Application Entity Titles may be site specific. Make sure that you check with the Filming Device Service Representative and the hospital network administrator to ensure you are using the correct AE Title for the destination filming device.**

9.2 Declaring DICOM Filming Devices on the CT System



NOTICE
Potential For
Data Loss

Empty all filming queues before modifying camera parameters.

This section contains procedures for recording important Camera setup information. Use the table(s) at the end of this section to record information from the setup screens.



- 1.) Click on the SERVICE DESKTOP icon.
- 2.) On the Desktop Toolbar select CONFIGURATION -> INSTALL CAMERA. The Install Camera window appears, along with a warning message pop-up box, to remind you that all filming queues must be empty before you begin to update or delete a camera.
- 3.) The Graphical User Interface displayed shows a list of cameras installed (See Figure 7-10). First, you must click OK in the warning message box. See Figure 7-9.

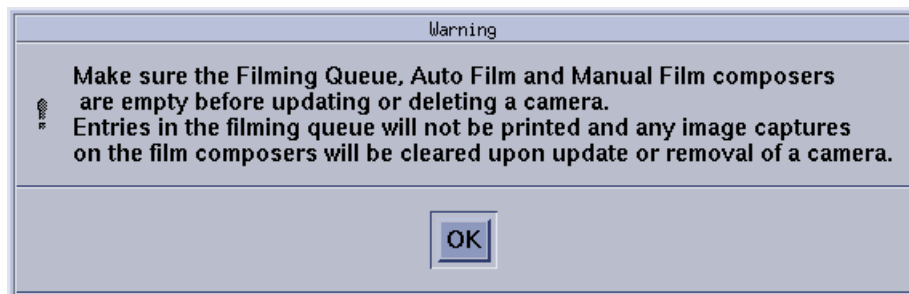


Figure 7-9 Warning Screen

- 4.) Now you are asked a series of questions.
- a.) To add a new camera, click the ADD button (See Figure 7-10).

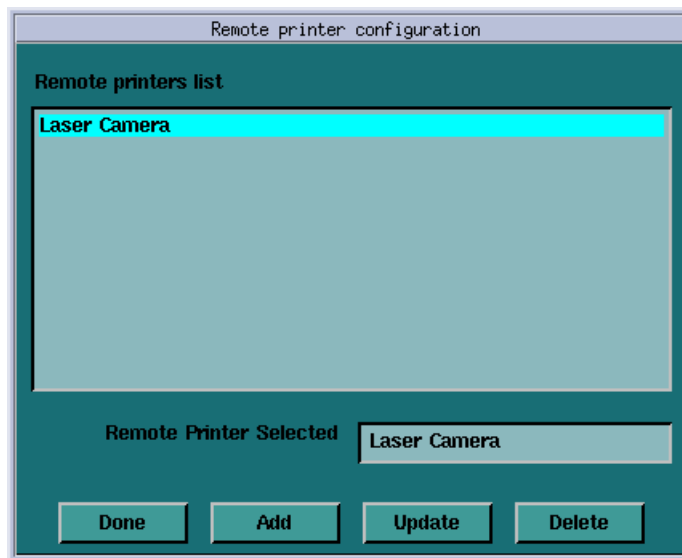


Figure 7-10 Printer Configuration GUI

- b.) A dialog window for the camera type (DASM/DICOM) appears. If no DASM is detected during the OC boot, the DASM button will be disabled (Figure 7-11). If a DASM is present and has not been detected, reboot the OC and run the camera configuration tool again.

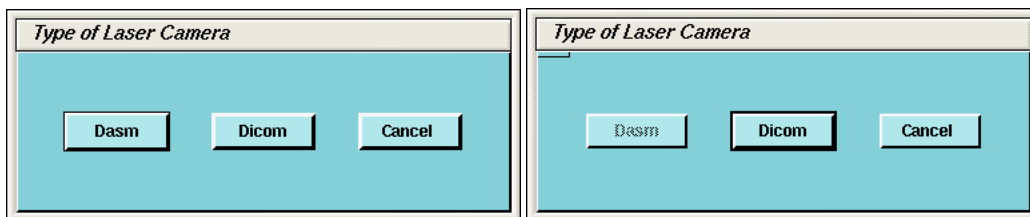


Figure 7-11 Dialog Box for camera Type

- 5.) To add a new laser camera, click DASM in the camera type dialog box. This brings up a list of available camera models. Select the appropriate model from the list and click SELECT (See Figure 7-12). Now you must configure it.

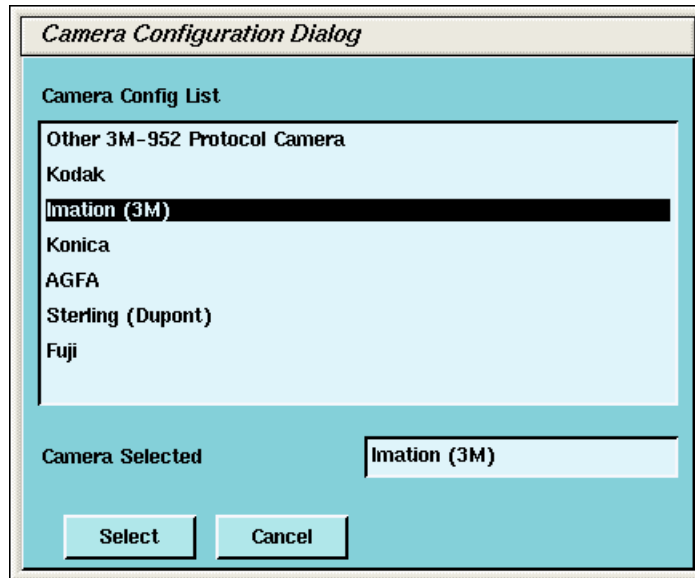


Figure 7-12 Camera Model Dialog (with DASM)

- a.) DASM Interface is automatically detected as either Analogue or Digital
- b.) Two Laser Options are available for laser cameras: SLIDES and ZOOM. Set this option only if the camera being installed supports slides and zoom. Setting the option allows it to be enabled or disabled at the application level.
- c.) Camera manufactures provide two (2) Magnification Type options for cameras. The SMOOTH resolution blurs the image, while the SHARP resolution makes the image pixels more pronounced. The default is smooth.

Comment:

To film good images, and have them look like images filmed by other GE CT products, the following camera settings are suggested:

Kodak: SMOOTH

Dupont/Sterling: SMOOTH

3M/Imation (Laser Camera): SHARP

3M/Imation (Dry View): SMOOTH

Agfa: SMOOTH

- d.) Select the appropriate File Format. Select ON from the drop down list boxes on the menu. Valid film formats are determined by the camera manufacture. IMATION for example, doesn't support 4x4, 2x4 or 1x2 and AGFA does not support 2x4) The DICOM print convention designates film formats by column and row (e.g. 12 on 1 film is 3x4).

When finished setting parameters, click on DONE and proceed to step 8.

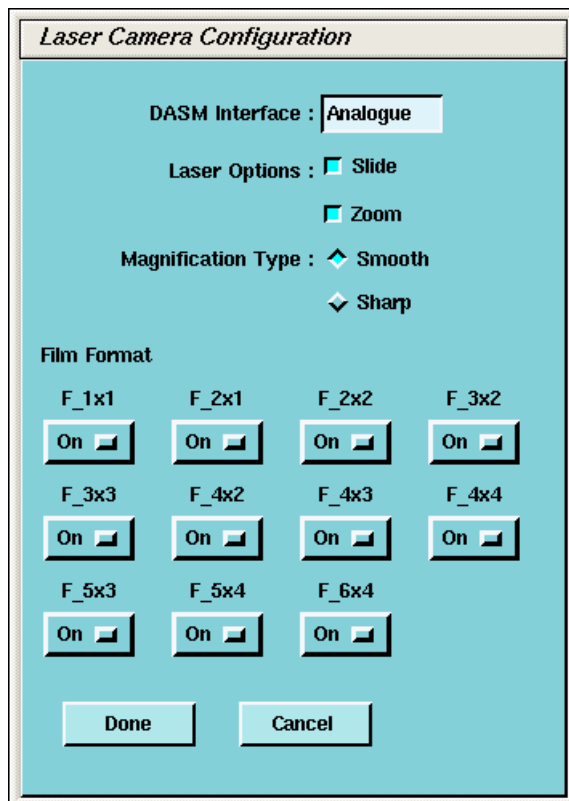


Figure 7-13 Laser Camera Configuration

- 6.) To add a new DICOM camera, click on ADD and then DICOM in the dialog box that appears.
 - a.) A list of camera models appears (See Figure 7-14). Select the appropriate model from the list and click SELECT. Clicking SELECT presets all the parameters to that models except the Network parameters.

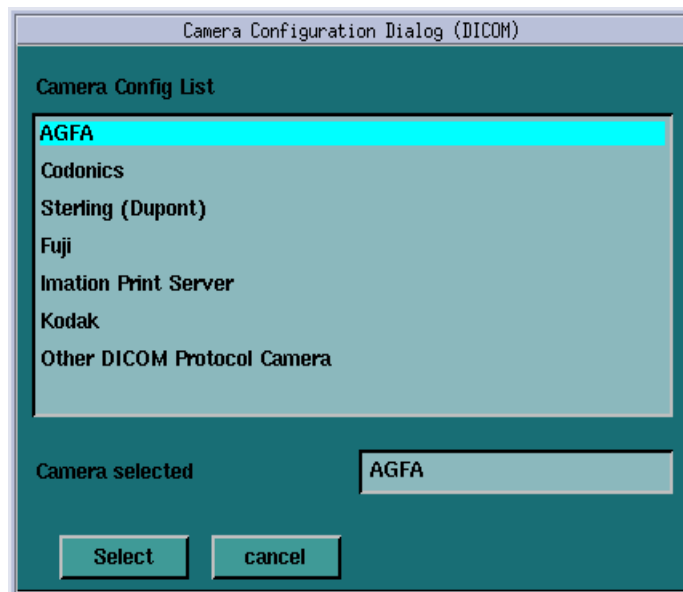


Figure 7-14 Camera Model Dialog (DICOM)

*Comment:
It's advised to
recheck the
preset
information with
the camera
manufacturer's
representative.*

- Selection of a different camera model clears the Image Quality parameters, because these are camera manufacture dependent.
- b.) Enter the Network Parameters (See Figure 7-15)
- >Device NameA unique name used to identify the camera.
 - >Host NameDICOM print server host name, as defined by the hospital.
 - >IP AddressDICOM print server IP address, as defined by the hospital.
 - >AE TitleDICOM print server application entity title, as defined by the print server. *You should consult the manufacturers DICOM Conformance Statement.*

Note: The Application Entity Title for the Camera may be site specific. Make sure that you check with the Camera Manufacturer's Representative and the hospital network administrator to ensure you are using the correct AE Title for the destination DICOM Print Camera.

- >TCP/IP Listen PortDICOM print server TCP/IP listen port, as defined by the server. *You should consult the manufacturers DICOM Conformance Statement.*
- >CommentsOptional comments used by the DICOM print server.

DICOM printer configuration parameters

Device Name:

Host Name:

IP address:

Application title:

TCP listen port:

Comments:

Medium Type:

Destination:

Film orientation:

Magnification type:

Standard film formats:

F_1x1 On	F_2x1 On	F_2x2 On	F_3x2 On
F_3x3 On	F_4x2 Off	F_4x3 On	F_4x4 On
F_5x3 On	F_5x4 On	F_6x4 On	

Figure 7-15 DICOM Camera Configuration

- c.) Medium Type specifies the type of film being used. Currently, only BLUE FILM and

CLEAR FILM are supported.

- d.) Set Destination to the final location for film output: either MAGAZINE or PROCESSOR.
- e.) Orientation selects film orientation. Only PORTRAIT is currently supported.
- f.) Set the Magnification Type. This parameter selects the algorithm used to interpolate pixels for proper film resolution. Set this parameter after consulting the camera manufacture to ensure optimal image quality. Choices are describe below:
 - >NoneNo interpolation. This option is not suported by all cameras.
 - >ReplicateAdjacent pixels are interpolated. This can result in "pixelized" images. *This algorithm is not normal preferred.*
 - >BilinearA 1st order interpolation of pixels. Results in images usuaully described as blurred. *This algorithm is not normal preferred.*
 - >CubicA 3rd order interpolation. Used with a large number of possible formulations. Camera manufactures define parameters called "smoothing type" to set coefficients used in this algorithm. The implementation of these "smoothing soefficients" is camera dependent.
- g.) Select the Standard Film Formats. Select the film format by choosing ON in the pull-down menu box located below each selection. See Figure 7-15. Valid film formats are set by the camera manufacture. IMATION for example, doesn't support 4x4, 2x4 or 1x2 and AGFA does not support 2x4) The DICOM print convention designates film formats by column and row (e.g. 12 on 1 film is 3x4).
- h.) After the camera has been configured, click the ADVANCED button. This creates the camera device file for you automatically and pops up the Advanced Parameters screen. See Figure 7-16.

Comment: For most Camera Manufacturers, the preferred selection is CUBIC.

The screenshot shows a dialog box titled "Camera Attribute Dialog". It contains several sections of controls:

- Image Quality Parameters** (Section Header)
- Configuration**: A text field containing "PERCEPTION_LUT=KANAMC".
- Maximum Density**: A text field containing "320".
- Minimum Density**: A text field containing "32".
- Smoothing Type**: A dropdown menu with "On" selected.
- Empty Density**: A dropdown menu with "Black" selected.
- Border Density**: A dropdown menu with "Black" selected.
- Film Size**: A dropdown menu with "<empty>" selected.
- Trim**: A dropdown menu with "Yes" selected.
- Priority**: A dropdown menu with "High" selected.
- Buttons**: "Done" and "cancel" buttons at the bottom.

Figure 7-16 Advanced Parameters (Camera Attribute Dialog)

- 7.) Advanced camera parameters control the image quality of films.

Note: For more information on the proper settings for these parameters, please see the Camera's DICOM Print Device Conformance Statement or the Camera Manufacturer Representative.

You may need to refer to a copy of the Conformance Statement as you are working with the Camera Manufacturer's Representative, to correctly set up the DICOM Print Camera I/Q and Time-out Settings.

- a.) **Configuration** - This parameter is camera manufacturer dependent as is typically used to specify the image contrast. The Configuration field may be up to 1024 characters long. The field will scroll automatically as text is entered. To review your entry, simply click and hold the middle mouse button, while the cursor is in the field, and drag the mouse towards the right (or left) as needed.

Note: Recommended Configuration Setting Values:

Agfa Drystar (MG3000)**PERCEPTION_LUT=KANAMORI (100)**
Imation Dryview (8700)**LUT=0,7**
Kodak Laser Printer 190**CS434\CN0\PD1.20**

- b.) **Smoothing Type** - Set Smoothing Type to ON, and input the selected value. This parameter is used when the magnification type is CUBIC. It represents the coefficient for the image resolution algorithm. This parameter is camera manufacturer dependent, and should be re-verified with your radiology department.

Note: **Recommended Smoothing Type Starting Values and Ranges:**

Agfa DryStar (MG3000)Start Value:**140**Range:**137 - 150**
Imation Dryview (8700)Start Value:**3**Range:**3 - 13**
Kodak Laser Printer 190Start Value:**Enhanced**Range:**Normal**

- c.) **Minimum and Maximum Density** - Used to set brightness of the images on film. The range of values is 0-4095, although the valid range for a specific camera is manufacture dependent. For **Maximum Density**, input the correct value into the text box. For **Minimum Density**, set it to ON and input the correct value in the text box.

Note: **Recommended Minimum and Maximum Density Starting Values:**

Agfa Drystar (MG3000)Min.:**20** or **23**Max:**300**
Imation Dryview (8700)Min.:**(Blank)**Max:**300**
Kodak Laser Printer 190Min.:**20**Max:**300**

- d.) **Empty Image Density** - This parameter sets the density for empty film viewports. Typically, BLACK is used but WHITE is an available option. The minimum and maximum density values are used as the representation.
- e.) **Border Density** - This sets the density for the border used around the film viewports. Typically, BLACK is used but WHITE is an available option. The minimum and maximum density values are used as the representation.
- f.) **Film Size** - Allows the system to specify a particular film size, if selected.
- g.) **Trim** - YES produces a white (clear) box surrounding each image.
- h.) **Priority** - This sets the print priority.
- i.) If you have completed entry of advanced parameters, click DONE.

CAMERA DATA TABLES: TO LOCATE INSTALL CAMERA INFORMATION: CLICK ON THE SERVICE DESKTOP BUTTON. ON THE DESKTOP TOOLBAR SELECT UTILITIES -> INSTALL -> INSTALL CAMERA. THE INSTALL CAMERA WINDOW APPEARS. SELECT EACH OF THE CAMERAS THAT ARE INSTALLED FROM THE LIST OF INSTALLED CAMERAS, AND CLICK ON UPDATE TO VIEW THE CAMERA'S SETTINGS. RECORD THE VALUES USED TO SET UP EACH CAMERA IN THE TABLES THAT FOLLOW. EXTRA TABLES

ARE PROVIDED FOR MULTIPLE CAMERAS

- Note: You can determine this information by looking at the contents of the following files:
- For a DASM Camera: `/usr/g/ctuser/app-defaults/devices/camera.dev`
 - For a DICOM Print Camera: `/usr/g/ctuser/app-defaults/devices/name.cfg` where, **name.cfg** is the camera device name from the printer configuration GUI.
- Example: `more <filename from above> ENTER`

DASM CAMERA #1 SETUP

GUI SETTING	SELECTIONS	VALUE
Camera Type	Model Type of Camera	
DASM Type	Digital or Analog	o Digitalo Analog
Options	Slides or Zoom	o Slideso Analog
Film	Smooth or Sharp	o Smootho Sharp
Film Format Available	1x1, 2x1, 2x2, 3x2, etc.	
Film Format Default	1x1, 2x1, 2x2, 3x2, etc.	

Table 7-1

DICOM CAMERA #1

GUI SETTING	SELECTIONS	VALUES
DICOM Camera Type	Model Type of Camera	
Film Format Available	1x1, 2x1, 2x2, 3x2, etc.	
Network Parameters	Host Name	
	IP Address	
	AE Title	
	TCP Listen Port	
	Comments	
	Destination	o Magazineo Processor
Special Set Up	Orientation	o Portraito Landscape
	Medium Type	o Blue o Clear
	Magnification Type	o None o Replicate o Bilinearo Cubic
*Advanced Parameters - IQ	Smoothing Type	o ON o OFF Value:
	Configuration	
	Minimum Density	o ON o OFF Value:
	Maximum Density	
	Empty Density (Black/White)	o Black o White

*To view Advanced DICOM Camera Settings, you must click on ADVANCED.

Table 7-2 DICOM Camera #1

GUI SETTING	SELECTIONS	VALUES
	Border Density (Black/White)	☐ Black ☐ White
	TRIM	☐ YES ☐ NO
	Priority	☐ HI ☐ MED ☐ LOW
	Film Size	

NOTES

*To view Advanced DICOM Camera Settings, you must click on ADVANCED.

Table 7-2 DICOM Camera #1

DICOM CAMERA #2

GUI SETTING	SELECTIONS	VALUES
DICOM Camera Type	Model Type of Camera	
Film Format Available	1x1, 2x1, 2x2, 3x2, etc.	
Network Parameters	Host Name	
	IP Address	
	AE Title	
	TCP Listen Port	
	Comments	
	Destination	☐ Magazine ☐ Processor
Special Set Up	Orientation	☐ Portrait ☐ Landscape
	Medium Type	☐ Blue ☐ Clear
	Magnification Type	☐ None ☐ Replicate ☐ Bilinear ☐ Cubic
*Advanced Parameters - IQ	Smoothing Type	☐ ON ☐ OFF Value:
	Configuration	
	Minimum Density	☐ ON ☐ OFF Value:
	Maximum Density	

*To view Advanced DICOM Camera Settings, you must click on ADVANCED.

Table 7-3 DICOM Camera #2

GUI SETTING	SELECTIONS	VALUES
	Empty Density (Black/White)	☐ Black ☐ White
	Border Density (Black/White)	☐ Black ☐ White
	TRIM	☐ YES ☐ NO
	Priority	☐ HI ☐ MED ☐ LOW
	Film Size	

NOTES

*To view Advanced DICOM Camera Settings, you must click on ADVANCED.

Table 7-3 DICOM Camera #2

DICOM CAMERA #3

GUI SETTING	SELECTIONS	VALUES
DICOM Camera Type	Model Type of Camera	
Film Format Available	1x1, 2x1, 2x2, 3x2, etc.	
Network Parameters	Host Name	
	IP Address	
	AE Title	
	TCP Listen Port	
	Comments	
	Destination	☐ Magazineo Processor
Special Set Up	Orientation	☐ Portraito Landscape
	Medium Type	☐ Blue ☐ Clear
	Magnification Type	☐ None ☐ Replicate ☐ Bilinearo Cubic
*Advanced Parameters - IQ	Smoothing Type	☐ ON ☐ OFF Value:
	Configuration	
	Minimum Density	☐ ON ☐ OFF Value:
	Maximum Density	
	Empty Density (Black/White)	☐ Black ☐ White
	Border Density (Black/White)	☐ Black ☐ White
	TRIM	☐ YES ☐ NO
	Priority	☐ HI ☐ MEDo LOW
	Film Size	

NOTES

*To view Advanced DICOM Camera Settings, you must click on ADVANCED.

Table 7-4 DICOM Camera #3

9.3 Troubleshooting Tips

The following is a summary of troubleshooting information for DICOM print that was gathered during software testing and validation of the DICOM print feature.

There is also a significant amount of additional troubleshooting procedures, and the theory of the DICOM print feature in the System Service Manual. Should you have problems installing a DICOM print camera, first read the information in System Service Manual.

ERROR BRINGING UP THE CAMERA INSTALLATION/FILMING APPLICATION

- **Symptom:** After creating/modifying the DICOM print camera the install camera interface does not come up and the filming application indicates it cannot bring up the film composer.
Cause: The configuration field contains invalid information such as a backslash (\) as the final character in the field or brackets ({}).
- **Solution:** The camera.dev file in ~ctuser/app-defaults/devices must be manually edited to remove the offending characters in the set configuration line. Invalid characters include \{ }

CONFIGURATION INFORMATION FIELD

- **Symptom:** Cannot view the entire configuration field (> 25 characters)
- **Solution:** Hold down the middle mouse button and move the field contents

NEED TO SET DICOM PRINT ATTRIBUTES NOT SUPPORTED BY SOFTWARE

- **Symptom:** User wants the white border around each image box ON/OFF permanently for this system and it cannot be set as the default for the camera.
- **Solution:** Using your favorite editor, add the following line to the camera.dev file located in ~ctuser/app-defaults/devices after the DICOM print device has been otherwise configured.

For Trim Off - set	TRIM NO
For Trim On - set	TRIM YES
- **Symptom:** DICOM print camera supports multiple film sizes and the user only wants to print if the film size is correct for BrightSpeed Select (14x17). [Otherwise the camera will queue the films or return an error causing the queue to pause (based upon the DICOM print camera specifications).]
- **Solution:** Using your favorite editor, add the following line to the camera.dev file located in ~ctuser/app-defaults/devices after the DICOM print device has been otherwise configured.

To force a 14x17 film size - set	filmSize 14INX17IN
----------------------------------	--------------------

NEED TO PREVENT DICOM PRINT ATTRIBUTES FROM BEING SENT TO DICOM PRINT CAMERA

- **Symptom:** Some DICOM print attributes are optional, and may result in fatal errors. For example, the Fuji camera does not support the Empty Image Density parameter for the film box.
- **Solution:** Using your favorite editor, add the following line(s) to the camera.dev file located in ~ctuser/app-defaults/devices after the DICOM print device has been otherwise configured.

- To prevent sending the Smoothing Parameter set	FB_Smooth FALSE
--	-----------------

-
- | | |
|--|-----------------|
| - To prevent sending the Border Density set | FB_Border FALSE |
| - To prevent sending the Empty Image Density set | FB_EID FALSE |
| - To prevent sending the Minimum Density set | FB_MinD FALSE |
| - To prevent sending the Trim Parameter set | FB_Trim FALSE |

ERROR TRYING TO CONNECT TO THE DICOM PRINT CAMERA

- **Symptom:** DICOM print server can be reached (ping), but Application error indicates "Unable to start filming interface" and the help message talks about running the install.dasm (Association Error)
- **Solution:** The system is unable to complete the association. Check the AE Title and the Port number of the DICOM print server and correct them through the Install Camera procedure.

FILM COMPOSER ERROR NOT USABLE

- **Symptom:** Film composer error says "unrecognized status - code 0"
- **Solution:** Review the log file, the attention and status windows. These areas have the correct filming status (for example, **film jam** and **supply empty**).

DEBUGGING CONNECTION ISSUES DIFFICULT

- **Symptom:** The timeouts for the DICOM print are very long, which means one needs to wait a long time before you know the application is not working.
- **Solution:** The timeouts for the DICOM print were setup to ensure that the system would work regardless of whether the DICOM print camera was on a LAN or a WAN halfway around the world. The DICOM print timeouts for the association and DIMSE classes (for example, **N-GET**, **N-DELETE**) can be modified within the DICOM print camera installation. They can be reduced down to 90 seconds.

DICOM PRINT ERROR ON N-GET TIMEOUT CONFUSING

- **Symptom:** When the N-GET timeout goes off, the error message in the prslog file will be "Could not get printer status, invalid command sequence for N-GET".
- **Solution:** When the user sees the above error they may want to consider that the issue may be an inactivity timer on the N-GET DIMSE service.

DICOM PRINT CAMERA SLIDE SUPPORT

- **Symptom:** Current implementation of DICOM print does not allow selection of slide format.
- **Solution:** Feature not currently supported. Possibly in future releases.

CONFUSION ON FILM FORMAT NOTATION

- **Symptom:** GEMS Laser Camera and DICOM Print film format notations are opposite.
- **Solution:**
 - GEMS Laser Camera film format notation has always been row x col (for example, 12 on 1 = 4x3)
 - DICOM Print Standard film format notation is col x row (for example, 12 on 1 = 3x4)

Section 10.0 Network Connections

BROAD-BAND

Broad-band is considered the standard network connection for BrightSpeed Select. (A dial-up modem is optional.) Broad-band connections should use one of the following Category 5 patch cables:

CAT Num	GE Part Num	Length
K9000WB	2215028-10	20 m
K9000KP	2215028-5	10 m
K9000JR	2215028-4	5 m
K9000WA	2215028-9	3 m

The CT system is connected to the network through the Console.

- A patch cable (not to exceed 10 feet) should be provided by the customer, and it is used to connect the console to a wall box. (See Notes on [Figure 7-17](#))
- Some customer-site units may require cable duct-work or conduit to route connecting network cables to the workstation, camera and console.
- The run from the hospital switch to the CT wall outlet must not exceed 290 ft. (88m). Bandwidth performance is degraded when the length reaches 300 ft. (91m) or greater.
- For the optional modem: **Two phone lines should be provided by the facility.** One line is for use with a modem and must be an analog line. The second line is a voice only line.

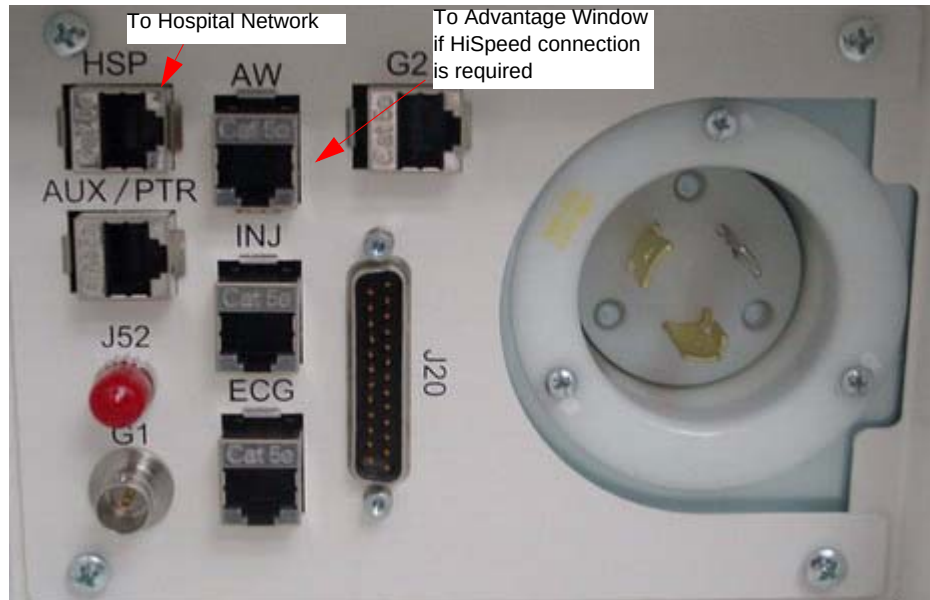


Figure 7-17 Console Rear Bulkhead

US PROCESS OVERVIEW

The United States network connectivity requirement for this product is broad-band. The US process relies on the Install Specialist to select a Customer Champion and identify an IT contact for the site. Together, those individuals then complete a site assessment to gauge what tasks are needed to fulfill the connection.

Anyone can contact the GE Connectivity team at 800.321.7937, Option #3, with questions.

CUSTOMER BROAD-BAND RESPONSIBILITIES

Provide GEMS Installation Specialist with an accurate site address, telephone number, contact name, and e-mail address for the:

- Customer Champion
 - Co-ordinate VPN activities between Radiology/Cardiology and the Information Technology (IT) departments
 - Act as a focal point in assuring site broad-band infrastructure meets GEMS requirements for connection as determined by a mutual assessment with the GEMS Connectivity team.
- IT Contact
 - Complete an equipment assessment with GEMS Connectivity team to determine site readiness for broad-band
 - Work with the Customer Champion to complete any identified infrastructure changes
 - Provide IP addresses for new CT equipment
 - Provide a VPN compatible appliance that will support the IPSec tunneling protocol and 3DES data encryption
 - To utilize an Internet Service Provider that supports static routing



GE Medical Systems

Remote Service Broadband - Customer Site Assessment

Site Name: _____ FE Name: _____
City, State: _____ FE Phone: _____
Date: _____ FE Email: _____

- | | Yes | No |
|---|--------------------------|--------------------------|
| 1. Does your site currently have a persistent (24x7) Internet connection? | ☐ | ☐ |
| 2. Is the GEMS Diagnostic Imaging equipment on the Local Area Network and will it be accessible to the Internet? | ☐ | ☐ |
| 3. Does your site have a VPN device today? | ☐ | ☐ |
| 4. Is the VPN device one of the models below? <i>If Yes, please select the model from the options below.</i> | ☐ | ☐ |
| ☐ a) Cisco Pix Firewalls ☐ g) Symantec (Raptor) firewalls
☐ b) Cisco Routers ☐ h) Firebox
☐ c) Cisco 3000 Series (Altiga acquisition) ☐ i) Linux S/WAN
☐ d) Checkpoint Firewalls Software
Version 4.1 and higher ☐ j) Sidewinder
☐ e) Nortel Contivity Software Version 3.2 or higher ☐ k) Netscreen
☐ f) Redcreek ☐ l) None
☐ m) Other _____ | | |

**If No, the GEMS Connectivity Support Team can help determine device compatibility.*

5. Does your VPN device support "triple DES" Encryption? ☐ Yes ☐ No
6. Has approval been given to install this VPN connection? ☐ Yes ☐ No
- Site Approver's Name _____
7. Provide your VPN Installer information, this is the person who will be contacted to schedule the VPN install.

Customer Installer Name: _____
Installer Telephone Number: _____
Installer e-mail address: _____

Notes: _____

Field Engineer needs to provide compatible system information:

System ID	IP Address	Gateway Address
_____	_____	_____
_____	_____	_____
_____	_____	_____
_____	_____	_____
_____	_____	_____
_____	_____	_____

[Additional System and IP Address Spaces Available on Page 2](#)

If you have questions or need assessment support contact your Zone Champ or:

Joe Gracz - HQ Support 1-262-524-5261
Joseph.Gracz@med.ge.com

Once you have completed both pages of this form, please send it to:

a) Judy Heyer judy.hey@med.ge.com
b) Judy Heyer Fax# 414-918-4707

[Use the send button on page 2](#)



GE Medical Systems

System ID List & Compatibility Matrix

MR	Software Version	CT	Software Version
Include all Signa(LX, Ovation, OpenSpeed and CV/I) MR's System Ids and IP Addresses in assessment.		HiSpeed CT/e	Minimum 6.0
		HiSpeed CT/e dual	Minimum 2.06
		HiSpeedNX/I	Minimum 5.5
		LightSpeed QX/I	Minimum 1.3
Prior to completing an InSite checkout systems must be upgraded to Version 9.0 or 9.1 unless shown in the compatible list below.		LightSpeed Plus	Minimum 2.1
Signa TwinSpeed	Release 9.0	LightSpeed Ultra	Minimum 3
Signa	Release 9.1	HiSpeed X/I	
Signa	Release 10.x	CT/I	Minimum 6.2
Signa CV/I	CNV4	Lxi, Dxi, Fxi, Zxi	Minimum 6.01
Signa 3T			
Signa Profile	Minimum 7.66		
Signa Contour	Minimum 7.66		
Magnet Monitor	Minimum 2.3 software		
Network Products	Software Version	Nuclear/Pet	Software Version
AW 4.0	Version 4.0 or Above	Advance	Minimum 5.1
AW 4.0P		eNTEGRA	Minimum 2.03

X-Ray

<u>Cardiac</u>	
INNOVA	All
<u>Mammography</u>	
Seno2000D	All
<u>Digital Radiology</u>	
Revolution XQ/i	Version 10.12.5 or above
Revolution XR/d	Version 18.0

Multi Vendor

MR	PSI	CT	PSI
Philips ACS 1.5T	MZP 400	Picker MX-TWIN	CZM 400
Philips NT 2000	MZP 401	Picker MX 8000	CZM 401
Philips NT 3000	MZP 402	Picker PQ	CZM 500
Philips NT ACS	MZP 403	Picker PQ 2000	CZM 501
Picker Outlook 0.23T	MZM 301	Picker PQ 5000	CZM 502
Picker Edge 1.5T	MZM 800	Picker PQ 6000	CZM 503
Picker Vista 1.0T	MZM 801	Picker PQS	CZM 504
Picker Eclipse 1.5T	MZM 900	Siemens AR STAR	CZS 100
Picker Polaris 1.0T	MZM 901	Siemens AR.C	CZS 101
Siemens Harmony 1.0T	MZS 100	Siemens AR.HP	CZS 102
Siemens Impact 1.0T	MZS 200	Siemens AR.SP	CZS 103
Siemens Impact Expert 1.0T	MZS 201	Siemens AR.T	CZS 104
Siemens SP 1.0T	MZS 400	Siemens PLUS-4 EXP	CZS 400
Siemens SP 1.5T	MZS 401	Siemens PLUS-4 POW	CZS 401
Siemens SP4000 1.5T	MZS 403	Siemens PLUS-4 VOL	CZS 402
Siemens Symphony 1.5T	MZS 500		
Siemens Vision 1.5T	MZS 600		

Field Engineer needs to provide compatible system information:


[illegible]

Once you have completed this form, please send it to:

Judy Heyer
Judy Heyer

email
Fax #

judy.heyer@med.ge.com
414-918-4707

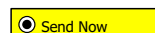

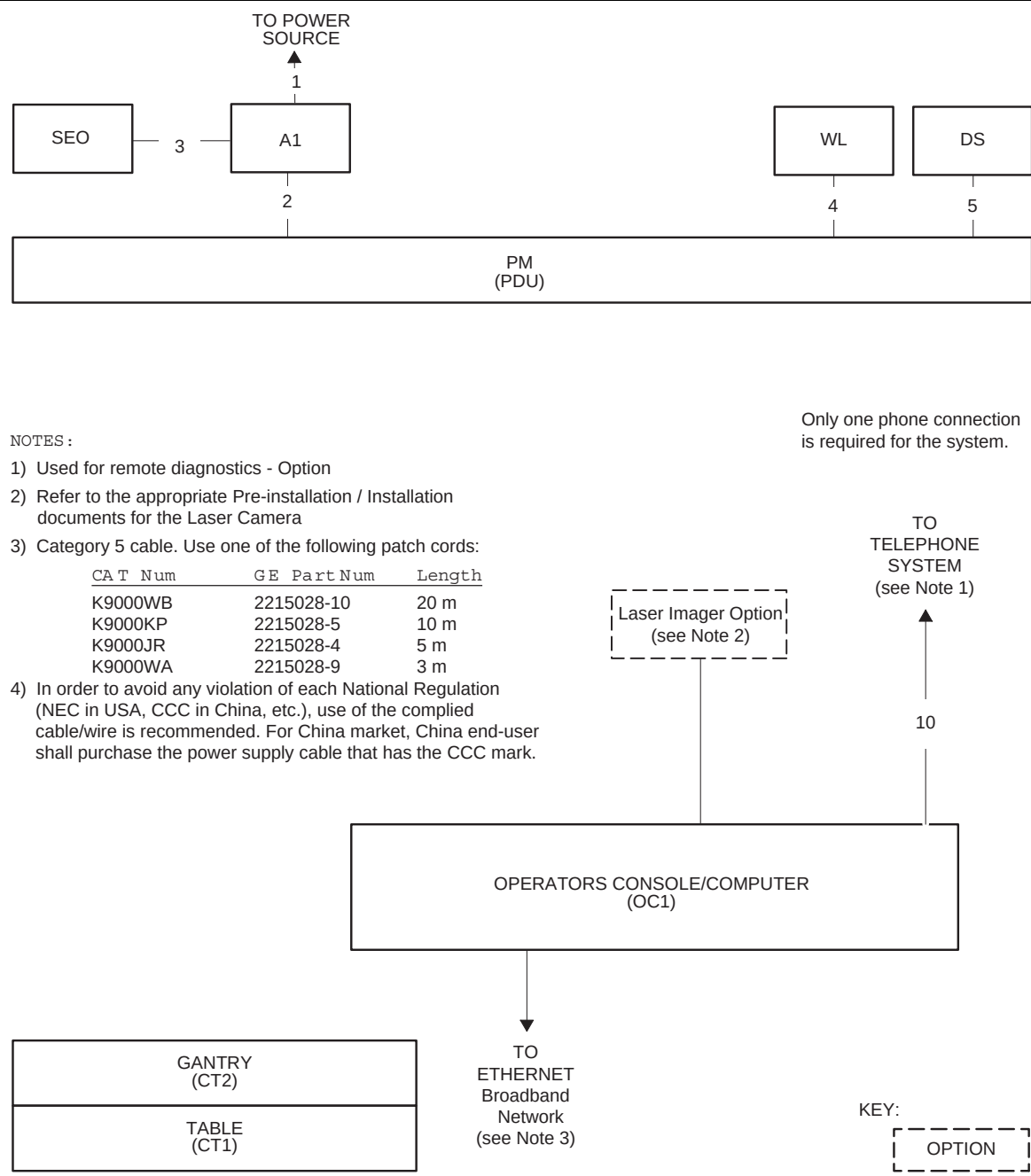


Figure 7-18 BrightSpeed Select System Interconnection Runs

Section 11.0

Modem Setup in IIP Configuration

If a phone Modem needs to be installed, follow below description to select correct modem type in "iipadmin config" GUI.

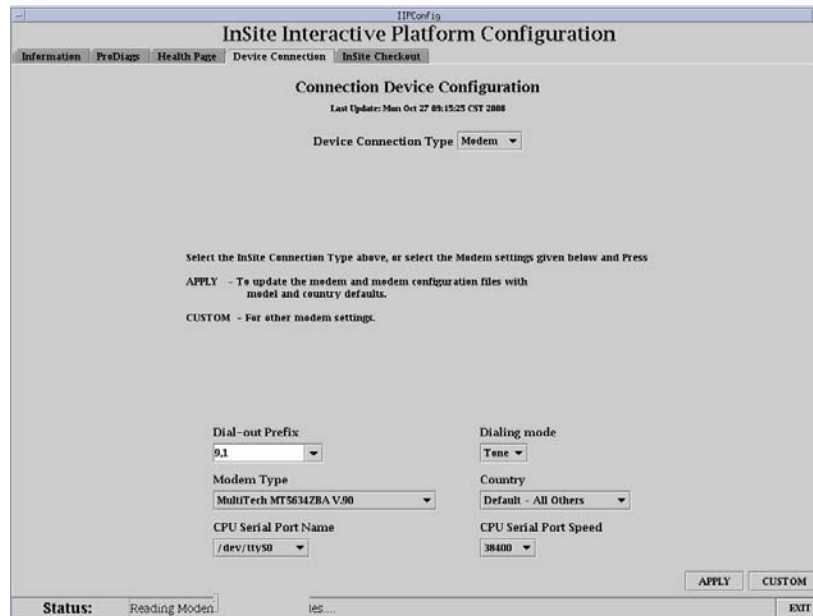


Figure 7-19 IIP Config GUI

For GOC4/LCGOC, a Series Port Modem will be used. At "CPU Series Port Name", Select /dev/ttyS0.

For All-In-One Console, see below:

If a USB Port Modem is used. at "CPU Series Port Name", Select /dev/ttyACM0.

If a Series Port Modem and USB Series converter cable is used. at "CPU Series Port Name", Select /dev/ttyUSB1.

Chapter 8

System-Level Safety Tests

Section 1.0

Patient Touch Leakage Test

1.1 Personnel Requirements

Required Persons	Preliminary Reqs	Procedure	Finalization
1 Engineer	10 mins	20 mins	N/A

1.2 Overview

This test should ensure that a patient undergoing a scan cannot touch or otherwise contact any conductive surfaces. Complete this procedure AFTER the installation of ALL options.

1.3 Preliminary Requirements

1.3.1 Tools and Test Equipment

Item	Qty
Standard service tool kit	1
Dale 600 / 601 Meter (from Tool Pool)*	1
Dale extended length leads*	1

*GE Healthcare performed validation of this procedure using the Dale 600/601 meter ONLY. Due to the unique nature of this meter, GE cannot guarantee the accuracy of this procedure if you use another meter.

1.3.2 Safety



WARNING



POTENTIAL FOR SHOCK

GROUND WIRES WILL HAVE GROUND CURRENT PRESENT WITH POWER “ON”. FOLLOW APPROPRIATE SAFETY PROCEDURES FOR WORKING WITH AN ENERGIZED SYSTEM.



NOTICE


Follow ALL required safety and PPE procedures customary for your organization, when working on this product.

1.3.3 Required Conditions

- Only trained service personnel should service the GE Scanner.
- GE Healthcare performed validation of this procedure using the Dale 600/601 meter ONLY. GE cannot guarantee the accuracy of this procedure if you use another meter.
- You must remove the Footswitch cover.

1.3.4 Test Conditions

- Test with the Table at maximum elevation and again with the Table at minimum elevation.
- Test to cover ALL points in an envelope described by Table travel from minimum to maximum extension, including the table extender.
- Test from both the head and foot of the Table -- e.g., test assuming that a patient may lie either HEAD-FIRST or FEET-FIRST.
- Test to ensure that all conductive surfaces are 6 feet (1.9 m) below and 7 feet (2.1 m) above the Table surface -- e.g., test assuming that a patient may lie either FACE-DOWN or FACE-UP and can touch or otherwise contact conductive surfaces both above and below the table.
- Test to ensure that a protective envelope, free of ANY and ALL conductive surfaces, extends at least 6 feet (1.9 m) around the patient. Surfaces of concern include:
 - I.V. poles and tray assemblies
 - Smart step monitors and stands
 - Table bearing rails, if accessible
 - The A1 Disconnect
 - The PDU
 - All other conductive surfaces

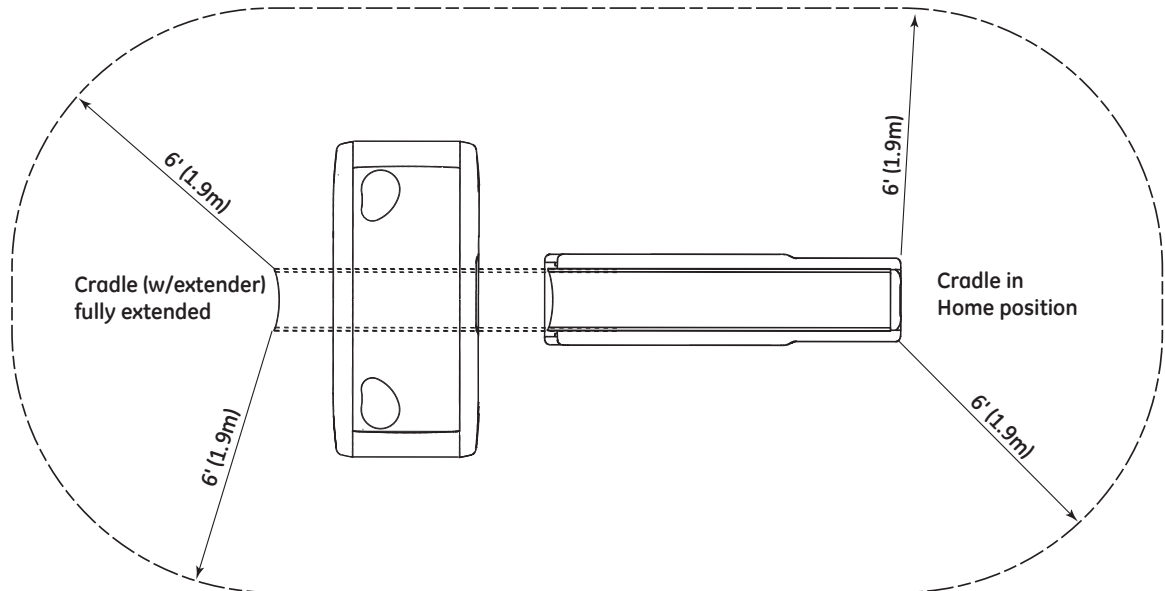

1.4 Procedure



NOTICE GE Healthcare performed validation of this procedure with the Dale 600/601 meter **ONLY**. Due to the unique nature of this meter, GE cannot guarantee the accuracy of these procedures when performed with any meter other than the Dale 600/601.

1.4.1 Cover Removal

- 1.) Move the Table to ISO elevation.
- 2.) Remove the Footswitch covers and the Gantry left side cover. Refer to [Appendix A: Removal & Installation of Covers](#) for additional information.

1.4.2 Leakage Current Testing

Note: Refer to the Dale 601 Operator's Manual for instructions on the use of the Dale 600/601 Meter for measuring leakage current or refer to [Figure 8-1](#) for a quick overview.

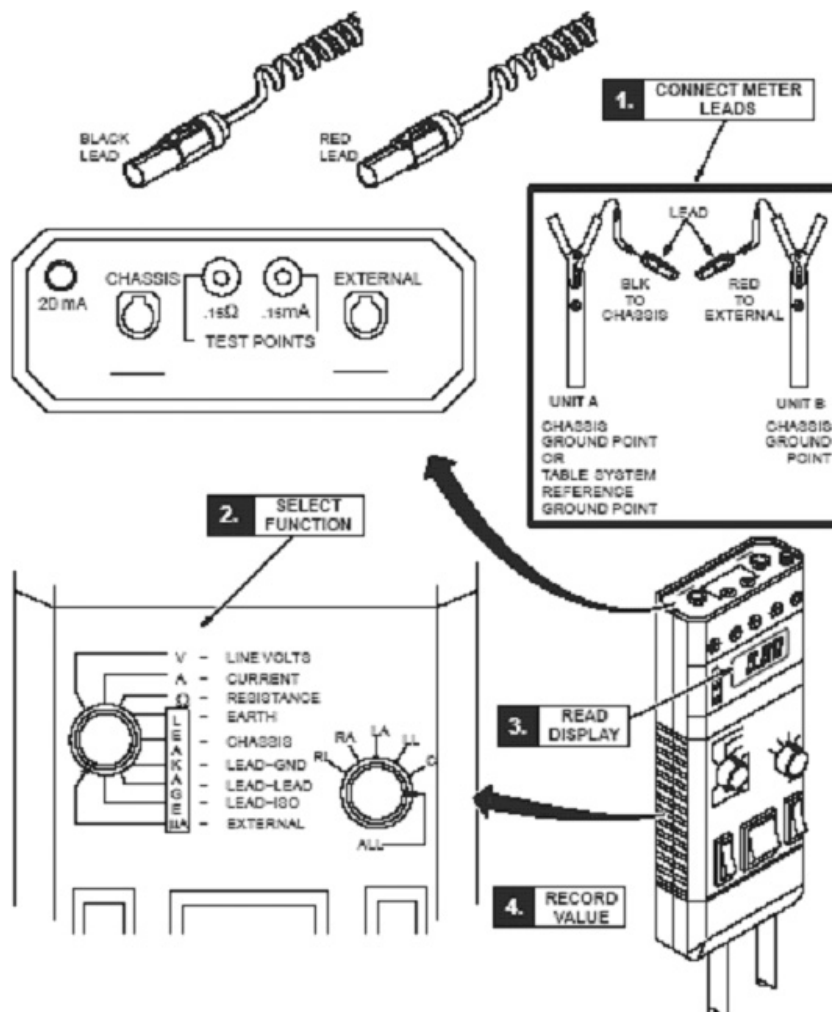


Figure 8-1 Using the Dale 600/601 Meter to Measure Leakage Current

- 1.) Plug the Dale 600/601 Meter into the outlet on the Gantry left side and confirm that the outlet is wired correctly by observing the three LED indicators on the meter.
- Note: Your meter must display a valid calibration sticker.
- 2.) Connect the meter leads to the meter as follows:
 - a.) Connect one end of the shorter black lead to the chassis plug and the other to the Table ground bus.
 - b.) Connect the red or longer black test lead to the external plug on top of the Dale 600/601 meter.
- Note: Some meters include a black lead and a red test lead, while other meters include two black leads. If your meter includes two black leads, the longer black lead is the test lead.
- 3.) Set the function switch on the Dale Meter to EXTERNAL and test that the meter is operational by touching the external lead to the meter's test terminal.
 - 4.) Keep the shorter black lead connected to the Table base ground bus and connect the red lead or longer black lead to the devices (components) under test, being sure to:
 - a.) Test with the power ON.
 - b.) Test between the system reference ground point (Table base) and the unit reference ground points (i.e. Gantry and Table, see [Figure 8-2](#)).
 - c.) Test all conductive surfaces and components within patient reach - within 6 feet (1.9 m) of the Table, and within 7 feet (2.1 m) above the surface of the Table, measuring at Table maximum travels. Refer to [Figure 8-2](#) for additional guidance.
- Note: At some sites, sinks and wall outlet cover plates may create areas of concern.
- d.) Test ALL optional components, such as in-room monitors, injectors, overhead monitor suspensions, and Table options.



NOTICE Beware of static discharge from the scan window, keypads, display, touch pads, and other plastic surfaces.

- 5.) Record the results of your tests in [Table 8-1](#), making certain that no leakage current exceeds the following limits:
 - * Critical care areas (invasive) - 10µA
 - * **General care areas - 20µA (Use this test criterion.)**
 - * Not intended for patient area - 50µA
- 6.) Complete the appropriate section of GE form e4879 to confirm the successful completion of this test. For additional information about this form, see Chapter 4, [Section 8.0 on page 189](#)

Table 8-1

Gnd Bus to	Install <20µA
Any Table ground points within the 6' range	
Any Gantry ground points within the 6' range	
Injector assembly metal surface	
Boom-in-Room metal surface	
Monitors or metal surface	
Sink or metal surface	
Installed Table accessories	
(other)	
(other)	
(other)	
(other)	

- 7.) Re-install the Footswitch cover and the Gantry left side cover. Refer to [Appendix A: Removal & Installation of Covers](#) for additional information.

PROCEDURE HINTS:

- Items with abnormally high or low measurements could indicate improper wiring, or loose or poor connections due to corrosion, painted surfaces, etc.
- High leakage could indicate a wiring error such as a neutral connected to the ground.
- Fluctuating ground currents could indicate a short, poor connection, or facilities ground problem causing leakage currents from other areas of the facility to flow through our system grounds.
- Refer to the illustrations below for visual depictions of the required measurements.

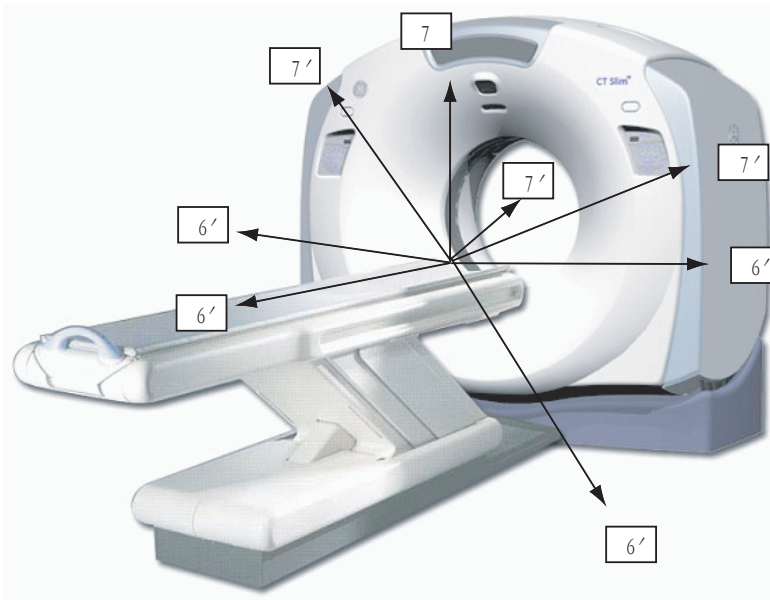


Figure 8-2 Measurement Points

1.5 Finalization

No finalization required.

Section 2.0

CT System Chassis Leakage Test

2.1 Personnel Requirements

Required Persons	Preliminary Reqs	Procedure	Finalization
1 Engineer	N/A	120 mins	N/A
1 Customer Electrician	N/A	120 mins	N/A

2.2 Overview

Unless required by your state, the CT System Chassis Leakage Test is an optional procedure.

2.2.1 Time

2 Hour labor on site, includes:

- Removal of covers
- Leakage test
- Re-installation of covers
- Recording of site data on the CT System Chassis Leakage Test Completion Form, if required by your state. Forward the completed form to your Project Manager of Installation.

2.3 Preliminary Requirements

2.3.1 Tools and Test Equipment

Item	Qty
Standard FE Service Tool Kit	1
GEH LOTO Kit	1
Dale 600 or 601 Leakage Current Meter	1

*GE Healthcare performed validation of this procedure using the Dale 600/601 meter ONLY. Due to the unique nature of this meter, GE cannot guarantee the accuracy of this procedure if you use another meter.

2.3.2 Safety



DANGER POTENTIAL FOR ELECTRIC SHOCK
GROUND WIRES WILL HAVE GROUND CURRENT PRESENT WITH POWER "ON".
THIS TEST IS PERFORMED WITH POWER "ON" - SYSTEM IN STAND-BY!



DANGER POTENTIAL FOR ELECTRIC SHOCK
SERVICING HARDWARE WITH POWER "ON."
EQUIPMENT SERVICE CAN ONLY BE PERFORMED SAFETY WITH THE MAIN
POWER "DISCONNECT" TAGGED AND LOCKED OUT.



WARNING Follow ALL required safety, PPE, and arc-flash procedures customary for your organization, when working on this product.

2.3.3 Required Documentation

Refer to Equipment Service section – LOTO PPE for detailed procedure.

2.3.4 Required Conditions

- Only trained service personnel should service the GE Scanner.
- Servicing hardware with power "On." Equipment service can only be performed safely with the main power "disconnect" Tagged- and Locked-Out.
- GE Healthcare performed validation of this procedure using the Dale 600/601 meter ONLY. GE cannot guarantee the accuracy of this procedure if you use another meter.
- You must remove the Footswitch cover and Gantry front, top, and side covers.

2.4 Procedure



NOTICE Follow LOTO and other safety procedures found in this manual before starting this procedure.

2.4.1 Cover Removal

- 1.) Remove the Gantry top and side covers. Guide rails on each side of the top cover will aid in removal. Refer to [Appendix A: Removal & Installation of Covers](#) for details.

Note: Fans located on the top cover make that cover heavier.

- 2.) Power down the console and follow the appropriate GEH Lockout/Tagout procedures.
- 3.) Remove the Gantry front cover, following the procedure in [Appendix A: Removal & Installation of Covers](#).

- 4.) Remove the Table footswitch top cover to gain access to the ground cables and ground bar.

Note: Do NOT disconnect ANY grounds at this time.

2.4.2 Electrical Connection Removal

- 1.) Confirm that all system grounds are securely attached to the system ground bus and NOT the Table base.
- 2.) Direct the electrician to remove all external electrical connections made during installation, including:
 - a.) Main system ground at PDU
 - b.) Power feeder flex connection at PDU
 - c.) Room door interlocks and room warning light connections
 - d.) Any and all other external ground connections to the system.
- 3.) Confirm the removal of all external gantry, table, console, and PDU connections, and have the electrician do the same.

Note: Some wires, such as the room warning light, may have external power and wire nuts, which should be installed to protect from arcing.

- 4.) Replace all covers removed in section 2.4.1 except for the Table footswitch cover.
- 5.) Follow the Lockout/Tagout procedure for re-energizing power and boot to application level. If not already done, remove the footswitch assembly top cover while the system boots.

2.4.3 System Ground Wire Testing

- 1.) Plug the Dale 600/601 Leakage Current Meter into one of the outlets on the gantry.
- 2.) Connect the meter leads to the meter as follows:
 - a.) Connect one end of the shorter black lead to the chassis plug and the other to the table ground bus.
 - b.) Connect the red or longer black test lead to the external plug on top of the Dale 600/601 meter.
- 3.) Set the function switch on the Dale 600/601 meter to EXTERNAL. Using the red or longer black external lead, touch the meter's test terminal to confirm that the meter is operational.

Note: For more information, refer to the Dale 600/601 Operator's Manual or see [Figure 8-1](#).

- 4.) With the system at application-level and all components functional, test the system ground wires as follows:
 - a.) Set-up the meter.
 - b.) Connect the shorter black lead to the system ground bus.
 - c.) Remove a system ground wire.
 - d.) Connect the red or longer black lead to the ground wire.
 - e.) Read the value on the meter.
 - f.) Replace the system ground wire.
 - g.) Record the results of your tests in [Table 8-2](#) for later use when filling-out the CT System Chassis Leakage Test Completion form.
 - h.) Repeat, testing all system ground wires ONE-AT-A-TIME. [Table 8-2](#) lists each system ground wire.

Note: The measured leakage current must not exceed 5 MA in any ground wire.

Table 8-2 System Ground Measurements

Components	Results	Date	Tester
PDU Leakage 1/0			
Gantry Leakage 1/0			
Console Leakage #2			
Table Leakage #10			
Option Grounds, if present			

2.4.4 System Power-Down

- 1.) After completing all tests, follow the Lockout/Tagout procedure to power down the system.
- 2.) Direct the electrician to re-install all electrical connections, conduits, cables, and wires removed in section 2.4.2 and to secure all connections per NEC code.
- 3.) Check to confirm that all connections are securely tightened.
- 4.) Re-install all removed system covers, except for the gantry right-side cover near the service switch panel.

2.4.5 System Power-Up and Test

Note: Be sure to complete this section of the procedure on-site.

- 1.) Follow the Lockout/Tagout procedure for re-energizing power.
- 2.) Turn-on the gantry service switches and power-up the console.
- 3.) Check that no cables remain in the gantry rotating path.
- 4.) Return cover dollies to storage areas.
- 5.) Check that the table controls and footswitch function properly.
- 6.) Re-test the system by completing a system functional scan. If installed, be sure to test the room warning light and the door interlock also at this time.
- 7.) Complete the CT System Chassis Leakage Test Completion Form, if required by your state, and forward the completed form to your Project Manager of Installation.

Note: Your product service information CD-ROM contains a copy of the CT System Chassis Leakage Test Completion form. For more information, see [GE and Regulatory Forms on page 189](#).

2.5 Finalization

No finalization required.

AppendixF

Additional Characterization Procedures

For information related to alignment, setup, and calibration procedures, please refer to the BrightSpeed Select System Service Manual, Direction 5181178-800.

AppendixG

System Configuration Data Sheets

REQUIREMENT

Record valuable system information in the data sheets that follow. Consult with your customer or network administrator to obtain the information. Understanding how the customer plans to use their CT scanner, and their network and filming expectation reduces the time required to reconfigure the system.

- [Table G-1 Manual Film Composer Options, on page 297](#)
- [Table G-2 System Network Configuration, on page 297](#)
- [Table G-3 Networking Application \(Image transfer\) Configuration, on page 298](#)
- [Table G-4 DASM Laser Camera Configuration, on page 298](#)
- [Table G-5 DICOM Print Camera Configuration, on page 299](#)
- [Table G-6 DICOM Print Camera Advanced Configuration, on page 299](#)

MANUAL FILM COMPOSER OPTIONS

MANUAL FILM COMPOSER OPTIONS

Slide Format (if available):

Greyscale:

Auto Printing:

Auto Clear Page:

Icon Labels:

Expose Order:

No. of Copies:

Table G-1 Manual Film Composer Options

SYSTEM NETWORK CONFIGURATION

SYSTEM NETWORK CONFIGURATION

	FIELD NAME:	SETENV NAME:	FIELD VALUE:
System Settings:	Service ID	SERVICE_ID	
	Hospital Name	HOSPITAL_NAME	
	Exam Number *	* ask customer or checklog	
	DAS Type	DASTYPE	
	PDU Type	PDUTYPE	
Network Settings:	Gateway Host Name	GATEWAY_HOSTNAME	
	Gateway IP	GATEWAY_IP	
	Gateway Net Mask	GATEWAY_NETMASK	

Table G-2 System Network Configuration

SYSTEM NETWORK CONFIGURATION			
	FIELD NAME:	SETENV NAME:	FIELD VALUE:
	Gateway Broadcast Mask	GATEWAY_BROADCAST	
	Suite Name	SUITEID	
Option	Network Printer IP Address		
Option	HIS Server IP Address		
Option	HIS Server AE Title		
Option	HIS server AE Port		
Option	CT Server AE Title		
Option	Connect Pro IP Address		

Table G-2 System Network Configuration (Continued)

NETWORK APPLICATION (IMAGE TRANSFER) CONFIGURATION

Record the network application (image transfer) configuration.

NETWORKING APPLICATION (IMAGE TRANSFER) CONFIGURATION				
AE TITLE OR HOST NAME	NETWORK ADDRESS	NETWORK PROTOCOL	PORT NUMBER	COMMENTS

Table G-3 Networking Application (Image transfer) Configuration

HOST ETHERNET ADDRESS

_____ : _____ : _____ : _____ : _____

CAMERA APPLICATION CONFIGURATION

Record the camera application configuration for the DASM or DICOM print camera.

DASM LASER CAMERA CONFIGURATION	
Camera Type:	

Table G-4 DASM Laser Camera Configuration

DASM LASER CAMERA CONFIGURATION	
DASM Type:	
Film Smooth/Sharp Setting:	
Options:	
Valid Film Formats:	
Default Film Formats:	

Table G-4 DASM Laser Camera Configuration (Continued)

DICOM PRINT CAMERA CONFIGURATION	
Camera Type:	
Host Name:	
IP Address:	
AE Title:	
TCP/IP Listen Port:	
Comments (Optional):	
Valid Film Formats:	
Default Film Formats:	
Destination:	
Orientation:	
Medium Type:	
Magnification Type:	

Table G-5 DICOM Print Camera Configuration

DICOM PRINT CAMERA ADVANCE CONFIGURATION	
Smoothing Type:	

Table G-6 DICOM Print Camera Advanced Configuration

DICOM PRINT CAMERA ADVANCE CONFIGURATION

Configuration:	
Minimum Density:	
Maximum Density:	
Empty Density:	
Border Density:	
Association Timeout:	
Session Timeout:	
N-Set Timeout:	
N-Action Timeout:	
N-Create Timeout:	
N-Delete Timeout:	
N-Get Timeout:	

Table G-6 DICOM Print Camera Advanced Configuration (Continued)

CONFIGURATION

Note: Type the text shown in **boldface**, and press the ENTER key on the keyboard.

NEXT STEPS

Resume installation following instructions in [Chapter 5: Restore System State \(Section 2.3\)](#).

AppendixH

Symbols

SYMBOL	PUBLICATION	DESCRIPTION
	417-5032	Alternating Current
	335-1	Three-phase Alternating Current
	335-1	Three-phase Alternating Current with neutral conductor
		Direct Current
	417-5019	Protective Earth (Ground)
	348	Attention, consult ACCOMPANYING DOCUMENTS
	417-5008	OFF (Power: disconnection from the mains)
	417-5007	ON (Power: connection to the mains)
		Warning, HIGH VOLTAGE
		Emergency Stop
		Type B

Table H-1 Symbols

SYMBOL	PUBLICATION	DESCRIPTION
	417-5339	X-ray Source Assembly Emitting
	417-5009	Standby
		Start
		Table Set
		Abort
		Intercom
		(on Operator Console) Power On: light on Standby: light off

Table H-1 Symbols

www.gehealthcare.com

