Nemoto Kyorindo

Dual Shot Alpha

Installation and Checkout Manual

INSTALLATION DOCUMENTATION

ADVANCED SERVICE DOCUMENTATION. PROPERTY OF NEMOTO
KYORINDO
USE OF THESE MATERIALS LIMITED TO AGENTS AND EMPLOYEES OF
NEMOTO KYORINDO OR OTHER PARTIES EXPRESSLY LICENSED BY
NEMOTO KYORINDO.
UNLICENSED USE IS STRICTLY PROHIBITED.

001-212134.00.03

Rev 3

© 2007, Nemoto Kyorindo Co., Ltd.
All Rights Reserved

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WARNING

(EN)

- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR, OR PATIENT FROM ELECTRIC SHOCK, OR FROM MECHANICAL OR OTHER HAZARDS.

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DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent write "Damage In Shipment" on ALL copies of the freight or express bill BEFORE delivery is accepted or "signed for" by a General Electric representative or hospital receiving agent. Whether noted or concealed, damage MUST be reported to the carrier immediately upon discovery, or in any event, within 14 days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this 14 day period.

To file a report:

- Call 1-800-548-3366 and use option 8.
- Fill out a report on <http://us44hdd21/sctg/InstallFulfill/InstallFulfillment.htm>
- Contact your local service coordinator for more information on this process.

CERTIFIED ELECTRICAL CONTRACTOR STATEMENT

All electrical Installations that are preliminary to positioning of the equipment at the site prepared for the equipment shall be performed by licensed electrical contractors. In addition, electrical feeds into the Power Distribution Unit shall be performed by licensed electrical contractors. Other connections between pieces of electrical equipment, calibrations and testing shall be performed by qualified GE Medical personnel. The products involved (and the accompanying electrical installations) are highly sophisticated, and special engineering competence is required. In performing all electrical work on these products, GE will use its own specially trained field engineers. All of GE's electrical work on these products will comply with the requirements of the applicable electrical codes.

The purchaser of GE equipment shall only utilize qualified personnel (i.e., GE's field engineers, personnel of third-party service companies with equivalent training, or licensed electricians) to perform electrical servicing on the equipment.

OMISSIONS & ERRORS

Customers, please contact your GE Sales or Service representatives.

GE personnel, please use the GE Healthcare iTrak Process to report all omissions, errors, and defects in this publication.

REVISION HISTORY

Rev	Date	Author	Reason For Change
1	March 5, 2008	J. Knipfer	Initial version.
2	April 2, 2008	J. Knipfer	Incorporate several corrections identified during the document review by GEHC. Add SSO and pressure test clarifications.
3	April 25, 2008	J. Knipfer	Converted document to Nemoto owned document per GEHC. Added Nemoto part number, revision history and cover picture.

TABLE OF CONTENTS

IMPORTANT PRECAUTIONS	1
REVISION HISTORY	9
TABLE OF CONTENTS	10
1- INTRODUCTION	11
1-1 Specifications	11
1-1-1 Dimensions	11
1-1-2 Weights	11
1-1-3 Electrical Ratings	11
1-1-4 Environment	11
1-1-5 Regulatory Approvals	12
2 - SYSTEM INSTALLATION	12
2-1 Unpacking	12
2-2 Assembly and Checkout Tools Required	13
2-3 Remote Stand Assembly	13
2-4 Console Assembly	14
2-5 Typical Room Configurations	14
2-6 Remote Stand Installation Planning	16
2-7 Ceiling Suspension Installation	17
2-8 Installation of Powerhead to Ceiling Mount	19
2-9 Overview of System Connections	20
3 - SYSTEM CONFIGURATION	21
3-1 User Settable Configurations	21
3-2 Service Mode Settable Configurations	21
4 - SYSTEM CHECKOUT	22
4-1 Basic Inspections	22
4-1-1 Power On	22
4-1-3 Powerhead Ram Speed Check	22
4-1-4 Lamp Indicator Check	22
4-1-5 Powerhead Rotation	22
4-1-6 Remote Stand Inspection	23
4-1-7 Syringe Clamp Inspections	23
4-1-8 Powerhead Motion Inspection	23
4-2 Performance Checks	23
4-2-1 Injection Time Check (Performance Injection #1)	23
4-2-2 Injection Volume Check (Performance Injection #2)	23
4-2-3 Injection Flow Rate Check (Performance Injection #3)	25
4-2-4 Injection Pressure Check (Performance Injection #4)	26
5 - INSTALLATION AND CHECKOUT COMPLETION	27
APPENDIX A – WARRANTY ACTIVATION AND INSTALLATION CHECKLIST	28

1- INTRODUCTION

This manual provides the instructions for installation and checkout of the Dual Shot Alpha contrast delivery system. It is intended and should only be used by those trained and qualified to complete the installation.

All local laws, regulations and requirements etc. shall be followed in conjunction with installation process defined herein.

1-1 Specifications

The following specifications are a subset of the complete product specification which can be located in the Operation Manual or Service and Parts Manual. This subset is provided as a quick reference during the installation process.

1-1-1 Dimensions

Powerhead	20.47 x 4.41 x 7.09 in.	520 x 112 x 180 mm
Control Console	4.21 x 10.63 x 9.76 in.	107 x 270 x 248 mm
Power Supply Unit	9.45 x 7.48 x 3.94 in.	240 x 190 x 100 mm
Remote Stand	25.43in.(dia) x 38.8 in.	646mm (dia) x 985.5 mm

1-1-2 Weights

Powerhead	15.43 lb.	9.2 kg
Control Console	8.16 lb.	3.7 kg
Power Supply Unit	13.23 lb.	6 kg
Remote Stand	37.7 lb.	17.1 kg

1-1-3 Electrical Ratings

Rated Voltage:	115 / 230 VAC + 10%
Frequency:	50 / 60 Hz
Power Consumption:	150 VA
Protection Type:	Class I
Protection Grade:	CF Type

FUSES

Power Supply Unit	5A 250VAC (power line 115 VAC) 2.5A 250VAC (power line 230 VAC)
-------------------	--

Syringe Heaters	F1 and F2, 0.5A rated, HM05 type
-----------------	----------------------------------

1-1-4 Environment

Transportation & Storage Temperature:	-40 to 70 degrees C
Transportation & Storage Humidity:	10 to 100% non-condensing
Transportation & Storage Atmosphere:	50 to 106 kPa

Operating Temperature:	10 to 40 degrees C
Operating Humidity:	30 to 75% non-condensing
Operating Atmosphere:	70 to 106 kPa

1-1-5 Regulatory Approvals

U.S. FDA 510k	K062168
cTUVus Mark (includes) UL 60601-1:2003 CAN/CSA-22.2 No. 601.1-M90	CU 72060872 0001
IEC 60601-1, -1-1, -1-4	120013030 003
CS Test Certificate 1	2013030 001

2 - SYSTEM INSTALLATION

2-1 Unpacking

The system is delivered in one package, unless the optional ceiling suspension has been ordered. Ceiling suspension systems are shipped in a separate individual package. Below is the bill of materials that will be included in the Dual Shot Alpha package. Follow the steps indicated when removing system components.

Table 1 – Bill of Materials

Remote Stand	Remote Stand	Ceiling Mount
Powerhead with remote stand arm	X	
Powerhead with ceiling mount arm		X
Control Console	X	X
Power Supply	X	X
Operation Manual	X	X
100mL Syringe adapter (2)	X	X
Control Console cable (3m)	X	X
Powerhead extension cable (20m)	X	X
Installation instructions	X	X
CF Card (preinstalled in Console)	X	X
CF Card storage box	X	X
Power cable	X	X
Cable hanger with screws	X	
Remote stand pole	X	
Remote stand base	X	
3mm Allen wrench	X	X
8mm Allen wrench	X	X
Anti-rotation screw		X
Warranty activation and Checkout form	X	X

2-2 Assembly and Checkout Tools Required

The following tools are required to assemble the Remote Stand assembly and connect the necessary cables.

ASSEMBLY TOOLS

- #2 Philips screwdriver
- #1 Philips screwdriver
- 3mm Allen wrench (included)
- 8mm Allen wrench (included with remote stand version only)


CHECKOUT TOOLS

- 10mL graduated cylinder
- 150mL graduated cylinder
- Pressure gauge
- Stopwatch

2-3 Remote Stand Assembly

Be sure to reference to Figure 2-1 and follow the assembly steps in order to assure proper assembly and operation of the Remote Stand assembly.

Figure 2-1

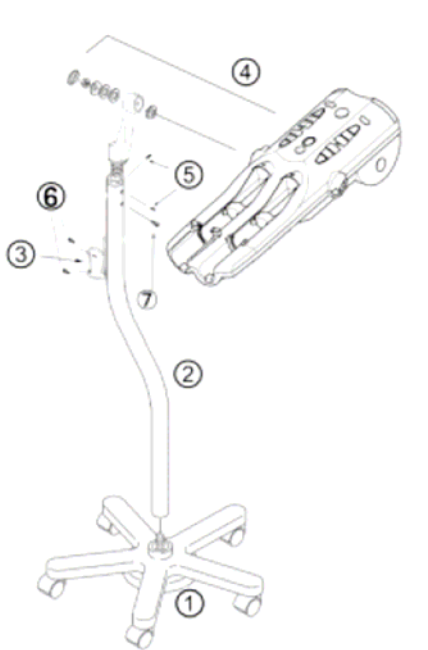Assembly steps:

1. Place the rolling base (1) on a level surface and insert the stand pole (2) into the opening. Turn the stand pole clockwise to thread onto the securing bolt until tightened.


2. Lay the stand pole and rolling base assembly on its side. Using the 8mm Allen wrench, tighten the hex bolt to secure the stand pole and rolling base. When completed, position assembly in the upright position.
3. Loosen the rotation limit screw (7), and locking screws (5) using 3mm Allen wrench. Insert the Powerhead arm into the top of the Remote Stand pole (2).
4. Tighten the locking screws (5) by inserting through the stand pole using the 3mm Allen wrench.
5. Tighten the rotation limit set screw (7).
6. Attach the Powerhead cable holder (3) using two M4 x 8mm screws (6) using Philips screwdriver.

2-4 Console Assembly

The Control Console is completely assembled before shipment. It is necessary however to remove shipping material before use. The CF Card which activates the Q-Select mode is installed in the Console before shipment and is protected from damage with packing material. After removing the Control Console from the shipping carton, remove the bubble rap material on the right-hand side of the console to expose the CF Card opening. Verify the CF Card is installed.

2-5 Typical Room Configurations

The following illustrations are provided to show typical CT room configurations and should be used for the purpose of the installation planning.

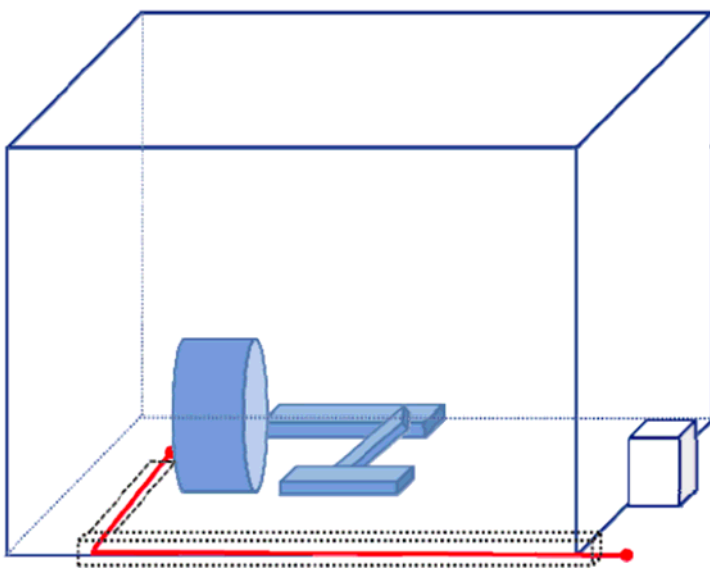This situation is ideal for a remote stand installation. A cable trough is provided in the floor of the CT suite. In this case the injector Powerhead cable may be routed through to the equipment room for connection to the Power Supply.

Be cautious when pulling cable through so as not to disrupt other cable connections.